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CERTIFICATE

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China Council for the Promotion of International Trade is China Chamber of International Commerce

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China Council for the Promotion of International Trade
China Chamber of International Commerce

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THIS IS TO CERTIFY THAT: the seal of SHENZHEN COSBEAUTY TECHNOLOGY CO., LTD. on the annexed DOCUMENT is genuine.

China Council for the Promotion
of International Trade

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Authorized

Signature:

Zhang Yuhong

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(Date: Aug. 27, 2021)

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Website for verifying the certificate: <http://www.rzccpit.com/validate.html>

INSTRUCTION FOR USE

Cosmetic device for photo
epilation BORK D680

Approved

Name: Lewis Huang

Position: CEO

Signature

Lewis Huang

NAME

Cosmetic device for photo epilation BORK D680

including:

1. Cosmetic device for photo epilation BORK D680- 1pc.
2. AC Adapter -1 pc.
3. Soft cloth for device care -1 pc.
4. Storage bag -1 pc.
5. Instruction for use with the warranty card -1 pc.

DESCRIPTION

Cosmetic device for photo epilation BORK D680 is used to remove hair on the skin surface.

DIRECTIONS FOR USE

BEFORE FIRST USE

- Read the user manual (operating instructions) carefully before using the device. Failure to follow the instructions for use specified in the manual may result in injury and damage to the user's property.
- Remove the device from the packaging.
- Remove and safely dispose of all packing materials.
- Use a razor to remove hair from the treated area. Do not replace this procedure with wax or

INSTRUCTION FOR USE

Cosmetic device for photo epilation BORK D680

epilatory cream.

Hair that is not completely removed from the treatment area can lead to burns, caking, and unpleasant odors and can be uncomfortable to perform. Burnt hair residues can adhere to the window or inner wall of the lamp operating window, adversely affecting performance and life.

After hair removal, the skin becomes sensitive and the user may experience discomfort during the procedures performed.

It is usually not necessary to remove hair before using the device if there have been 4 or 5 sessions of its use.

Also, it is not necessary to remove hair before use if it does not grow or its growth does not become evident.

- Cleanse the skin and make sure it is completely dry and free of visible hair and cosmetics containing oils (make-up, deodorants, perfumes, sunscreen and sunscreen).
- Check that there is no hair, dust or fibers on the appliance, especially on the working window.

Table 1. Table of hair color, skin tone and recommended intensity level

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PROCESSING LEVELS

If you are using the device for the first time, start at level 1. After adapting to the senses, you can set the device to a higher level.

The higher the level, the stronger the processing effect will be. It is important to pay attention to the state of the body during use. If you feel uncomfortable, then reduce the level of use.

The change in power is required to treat different body parts and skin tones. For example, the skin tone on the outside of the arm is darker than on the armpits because the outside is more exposed to the sun. A darker skin tone absorbs more rays of light, making it easier to irritate such skin. Therefore, it is recommended to set a lower intensity for dark skin treatments.

EXPLOITATION

INSTRUCTION FOR USE

Cosmetic device for photo epilation BORK D680

For first use, select a body part, such as an arm, and set the minimum treatment level.

After trial treatment, wait 24 hours and check for negative reactions.

This procedure will allow you to adapt to the light treatment procedures, and then you can choose the appropriate intensity according to your actual experience of using the device on your skin.

- To turn on the device, connect the device to the power adapter and press the power button. The device starts up and its screen turns on. The screen will display the processing levels, with the default set to level 1, auto shut off.
- To select epilation level, quickly press the power button to select a level between 1-5.
- Place the device at a 90 ° angle, place the working lamp window directly on the skin; the appliance will begin to emit light, which will stop when you move the appliance away from your skin.

Start moving the device, keeping the light source perpendicular to your skin, the device will continuously emit light.

- After use, press and hold the OFF button until the appliance is turned off completely.

RECOMMENDATIONS FOR USE

The device does not have a separate button for turning on the light source. If it detects skin, it will automatically start emitting light in the format of rapid successive flashes.

It is forbidden to place the working window of the lamp on the skin when starting the device, and also to touch the skin with the device immediately after starting it until the required treatment level is set.

Please do not move the device multiple times over the same skin area.

The instrument will turn off automatically three minutes after there is no activity.

The sound of the cooling fan can be quite noticeable. This is normal as this element is designed to reduce the temperature of components that become very hot during operation.

The sound emitted each time the light source blinks is normal, a higher intensity in the energy supply will increase this sound.

A slight tingling sensation and temporary redness are normal. They appear more obviously around hair growth areas.

PROCESSING PLAN

The total hair growth cycle is 18-24 months, during which continuous treatment with a photoepilator at specified intervals may be required to achieve a stable long-term result. The effect of the procedure varies depending on the condition of the body, treatment method,

INSTRUCTION FOR USE

Cosmetic device for photo epilation BORK D680

dy part, hair color and hair density.

Recommended treatment frequency during the full hair growth cycle:

The first four procedures should be performed once every two weeks.

The fifth-seventh procedure is performed once every four weeks.

- Treatment after the eighth procedure can be carried out regularly in accordance with the growth and condition of the hair. It is proposed to carry out the procedures once every four weeks or once every eight weeks in accordance with the state of health, while the break between the procedures must be at least two weeks.

If you want to treat another part of the body immediately after using the device in any area, then you need to check the cleanliness of the working window of the lamp.

For best results, it is recommended that you clean the device after treating each part of the body.

We recommend using a mirror when treating the underarms, upper lip and chin for maximum precision and comfort.

The lamp has a total life of 300,000 flashes. Use must be stopped when the flash counter reaches "000000".

SKIN CARE RECOMMENDATIONS

The treated skin should be moisturized with a special product or lotion. If the skin becomes red or itchy, cool it down with a damp towel.

Do not use cosmetics with fragrances and scrubbing ingredients on the treated skin area within 48 hours after the procedure.

Within 24 hours after hair removal, the skin area must be protected from UV rays and sunscreen must be applied during the treatment period.

If the skin remains red after 24 hours after treatment, it must be additionally protected from ultraviolet rays until the color returns to normal.

See your doctor immediately if you experience any problems.

CLEANING AND MAINTENANCE

- Switch off the device or unplug the power cord from the outlet.
- Let the appliance cool down and clean it with a damp soft cloth.

Clean the working lamp window, skin sensor and instrument surface. A dry cotton swab can be used to clean the edges and corners of the appliance. If the device cannot be cleaned effectively with a damp cloth, alcohol may be applied to the surface.

INSTRUCTION FOR USE

Cosmetic device for photo epilation BORK D680

To avoid dropping and damaging the unit, place it in the storage case after use.

OPERATING, STORAGE, TRANSPORTATION CONDITIONS AND SHELF LIFE

Ambient temperature	5 - 30°C
Relative humidity	Not more than 65%
Air pressure Range:	86 kPa – 106 kPa

It is recommended to store the device in a cool dry place, isolated from places where acids and alkalis are stored.

Protect the device from significant temperature fluctuations and direct sunlight.

The shelf life of the product from the date of sale is 5 years.

Service life - 3 years from the date of sale.

The warranty period is 1 year from the date of sale.

The lamp has a total life of 300,000 flashes. Use must be stopped when the flash counter reaches the

"000000".

TROUBLE-SHOOTING

Failure	Possible reason	Decision
Light source cartridge won't turn on	Window is not correctly positioned against the skin	Place the window on the skin correctly
The device does not turn on	The device is not connected correctly	Disconnect the instrument and then reconnect it to the power adapter
During use, there is a burning sensation and pain	The selected processing level is too high	Press the power button to select a lower level
	There is hair on the skin	Remove hair before use

INTENDED PURPOSE

Medical device MD «Cosmetic device for photo epilation BORK D 680» is intended to remove hair on the skin surface.

INSTRUCTION FOR USE

Cosmetic device for photo epilation BORK D680

INDICATIONS

to remove hair on the skin surface.

CONTRAINDICATIONS

It is forbidden to use the device by persons suffering from:

- tuberculosis,
- fever,
- hypertension,
- infectious diseases,
- skin diseases,
- having malignant tumors of various kinds,
- a state of acute allergic reaction,
- blood diseases, diseases of internal organs;
- persons suffering from atopic dermatitis;
- persons with sensitive skin,
- increased fatigue or poor physical condition;
- after drinking alcohol,
- after taking sleeping pills;
- women during pregnancy, breastfeeding, menstruation;
- persons suffering from hemorrhage syndrome,
- epilepsy,
- persons with pacemaker
- diseases of photosensitivity, for example, polymorphic photodermatosis, urticaria, and lupus erythematosus;
- persons suffering from angioedema caused by long-term use of steroid hormones;
- to persons suffering from hepatitis;
- persons with suppressed immune response (including HIV or AIDS);
- people receiving radiation therapy or chemotherapy within the last three months.

INSTRUCTION FOR USE

Cosmetic device for photo epilation BORK D680

THE EFFECTS

tingling sensation
redness
burns
blisters,
persistent irritation
allergic reaction

⚠ PRECAUTIONS

- Store the device at a sufficient distance from sources of flame.
- Keep the device out of the reach of children.
- Do not touch the device with sharp objects or objects with high temperature.
- Do not twist the power cord around the power plug.
- Do not use accessories such as power adapters that are not supplied with the machine.
- Do not use the appliance in a humid environment, such as a bathroom, to avoid electric shock.
- Do not use the appliance in the following cases: the power cord or plug is damaged, smoke comes out of the openings of the appliance, a burning smell is felt during use, the appliance does not work or is damaged, external elements and the filter are damaged.
- Do not allow water and foreign bodies to enter the device body.
- Do not use the appliance on artificially tanned skin as this may irritate the skin and cause discomfort.
- Do not use the device on the skin after sunburn within 24 hours; Avoid exposing treated skin to the sun for 24 hours as this can irritate and damage the skin.
- Keep the device away from your eyes to avoid damage.
- Do not treat areas of skin with tattoos, pigmentation (for example, freckles, moles or nevus). Using the device on skin with excessive melanin content (dark skin) may cause side reactions such as burns or redness.
- Disassembly, repair and maintenance of the device is possible only in an authorized service center.
- Do not drop the device to avoid damaging it.
- The device is intended for personal use, for sanitary reasons it is not recommended to use it by several persons.

INSTRUCTION FOR USE

Cosmetic device for photo epilation BORK D680

The device is suitable for use by people with different skin types, except for dark-skinned and very dark-skinned. It is also suitable for use with most hair types, except for blonde, gray and light blond hair. If your skin or hair color does not match the specified requirements, it is recommended to refrain from using this device.

For the first use of the device, select several areas of the skin on the arm for testing, monitor the condition of the skin for 24 hours. Stop using the appliance in case of discomfort or allergic reaction.

- Avoid multiple use of the device on the same skin area, as this increases the likelihood of negative side effects. If you experience skin reactions with the appliance, such as burns, blisters, persistent irritation, stop using the appliance immediately and consult a doctor.
- Do not use the device on the following areas of the body: around the eyes, eyelids, inner ear, nostrils, stab wounds, areas that have undergone plastic surgery or other interventions, oral cavity and mucous membranes, areas with eczema, inflammation from excessive burns, acne and its varieties, with inflammation resulting from suppuration of the skin, itching fever; moles, bulging pigmented nevi and freckles, areas of large vessels, areas with a darker skin tone, with scars or skin abnormalities, large veins. It is forbidden to handle the Adam's apple, nipples, areoles of the mammary gland, navels, genitals, anus and surrounding areas, warts, tattooed areas, areas with varicose veins or dilated blood vessels, areas of contact with silicone implants, stimulants, subcutaneous injection holes (for an insulin syringe), piercing sites and so on.
- Do not use the device after taking any of the following medications: skin has been or has been treated with AHA and BHA acids, or has been partially treated with Roaccutane and azelaic acid within the last 30 days;
- Acutane® A or Roaccutane® A has been used orally or topically within the past six months;
- the user has taken or is taking medications that make the skin particularly sensitive to damage, injury or exposure to irritants; • the user has taken or is taking pain relievers (as these usually dull the feeling of warmth);
- the user has taken a photosensitizer or a similar medication (refer to the instructions enclosed in the medication box, do not use the device if there is a risk of photosensitization or phototoxicity reactions, or if the instructions prohibit exposure to the sun after taking medication);
- the user has taken or is taking anti-coagulant (blood thinning) drugs, including large doses of aspirin (drugs must be stopped at least one week before the session to facilitate the elimination of drugs from the body);
- the user has taken or is taking immunosuppressants.

LIGHT EXPOSURE WARNINGS

It is unacceptable to direct the light source into the eyes.

WARNING! Do not use the hair remover near your eyes

INSTRUCTION FOR USE

Cosmetic device for photo epilation BORK D680

on the forehead.

eye safety

Do not look directly into the lamp window or attempt to activate the device toward the eyes. To protect the safety of your eyes and those of individuals nearby, the device will only activate when skin sensor is placed against the skin. If at any time the device activates when lamp window is NOT in contact with the skin, STOP use immediately and contact your local representative. Possible eye injury (potentially leading to loss of vision) or skin injury if instructions are not followed. Protect the eye from exposure. Do not use the appliance over the eye lids or close to the eye.

REPAIR AND TECHNICAL SERVICE

Disassembly and maintenance of the device is possible only at an authorized service center.

TROUBLESHOOTING MEASURES

- If a malfunction is found in the operation of the device, immediately stop using it and take it to the nearest BORK service center for testing, replacement or repair.
- For safety reasons, all repair and maintenance work on the device should only be performed by representatives of the authorized service centers of the company.

Service center address

OOO "Premium Service", Moscow, Bolshaya Sadovaya st., 10.

Up-to-date information on the availability of service centers in your region on the website bork.ru or by phone 8 800 500 88 99.

CLASSIFICATION

Electric shock protection	Class II
Working parts class	Type BF applied parts (according to IEC 60601-1)
Protection against ingress of liquids:	IP22

GUARANTEED BY MANUFACTURER VALUES OF MAIN PARAMETERS, CHARACTERISTICS (FEATURES) OF MEDICAL DEVICE

Where applicable, a deviation of $\pm 10\%$ from the nominal characteristics is possible, unless otherwise indicated.

INSTRUCTION FOR USE

Cosmetic device for photo epilation BORK D680

Main Technical Parameters of the medical device guaranteed by manufacturer

Characteristics	Value
Product Name	Cosmetic device for photo epilation BORK D 680
Model No.	D 680
Dimensions	197.50*70.0*43.0mm $\pm 10\%$
Main Unit Weight	240,5 g $\pm 10\%$
Treatment Area	3.12cm ²
The power of the light flow	2,8 J/cm ² - 4,1 J/cm ²
Spectrum	510nm
Energy Level	5 levels (from 8 till 15 J) 1 level : 9.5J ± 1.5 J 2 level : 10.5J ± 1.5 J 3 level : 11.5J ± 1.5 J 4 level : 12.4J ± 1.5 J 5 level : 13.1 ± 1.5 J Correspondence between Energy level and skin and hair type see table 6
Power Supply	AC 100~240V,
Rated frequency	50/60Hz
Output Voltage/Current	DC 12V/4A
Rated Power	48W
Lamp	Lamp type – quartz Xenon 30 W Manufacturer: Shenzhen Senler Photon Technology Co., Ltd (China) lamp replacement is not available

The device is operated on the basis of software:

firmware for the device

INSTRUCTION FOR USE

Cosmetic device for photo epilation BORK D680

Software version: V1.0.11

WARRANTY OBLIGATIONS

The company is under warranty for one year from the date of sale of this product. The service life established by the company for this product is 3 years from the date of sale, provided that the product is used in accordance with the rules and recommendations set forth in the instructions for use attached to it and the applicable technical standards.

In order to avoid misunderstandings, we kindly ask you to carefully study the instructions for use, conditions of warranty obligations, check the correctness of filling out the warranty card. The warranty card is valid if it contains correctly and legibly indicated: model, serial number of the product, date of sale, clear seals of the seller, signature of the buyer. The model and serial number on the product must match those indicated in the warranty card.

In case of violation of these conditions and in the event that the data specified in the warranty card is changed, erased or overwritten, the warranty card is invalidated. If the date of sale cannot be established, in accordance with the law "On Protection of Consumer Rights", the warranty period and service life are calculated from the date of manufacture of the product.

The date of manufacture is encoded in the serial number of the product,

for example: 1st and 2nd signs - week of manufacture; 3rd and 4th signs - year of manufacture; the last 5 digits are the serial number.

WARRANTY CARD:

Medical device	Cosmetic device for photo epilation BORK D 680
Model	BORK D680
Serial number	_____
Date of purchase	_____
Information about the seller (name, legal address)	_____ _____

INSTRUCTION FOR USE

Cosmetic device for photo epilation BORK D680

TERMS OF SERVICE

The manufacturer reserves the right to refuse to satisfy consumer requirements for warranty obligations in case of non-observance of the operating conditions, installation of the product, deliberate or careless actions of the buyer (consumer) or third parties that caused damage to the product.

The warranty does not cover:

• elements and parts with natural wear. At the same time, natural wear is understood as the consequences of the operation of the product, which caused a deterioration in its technical condition and appearance due to prolonged use of this product.

The warranty does not apply to products damaged by:

- misuse, neglect, violation of storage rules, incorrect connection, non-observance of the attached user manual (operating instructions);
- ingress of foreign objects, substances, liquids, insects or animals inside the product body;
- exposure to high and low external temperatures on non-heat-resistant parts of the product;
- mechanical influences;
- unskilled repairs and other interventions that entailed changes in the design of the product;
- use of the product for industrial and / or commercial purposes;
- damage to electrical cords;
- use of cleaning agents not specified in this user manual (operating instructions);
- the use of non-standard (non-original) consumables, accessories, spare parts;
- connection to a network with a voltage different from that specified in the user manual (operating instructions) of the product, or with a mismatch between the parameters of the supply network and state standards of the Russian Federation;
- force majeure circumstances (fire, flood, lightning).

Replacement of faulty parts (parts, assemblies, assembly units) in the product during the warranty period does not lead to the establishment of a new warranty period for the entire product or for replaced parts.

MATERIALS

Body – ABS

POTENTIAL USERS

User/Operator:

INSTRUCTION FOR USE

Cosmetic device for photo epilation BORK D680

adult (the operator, has received education and intensive reading ability, who can read instructions);

Health: Good or (patients can operate themselves under the doctor's supervision).

Nationality or Race: Unlimited.


Able to identify body parts.

Place of application: at home or in a hospital.

LABELING

1. Rating plate (identification label)

Located on the back of the device

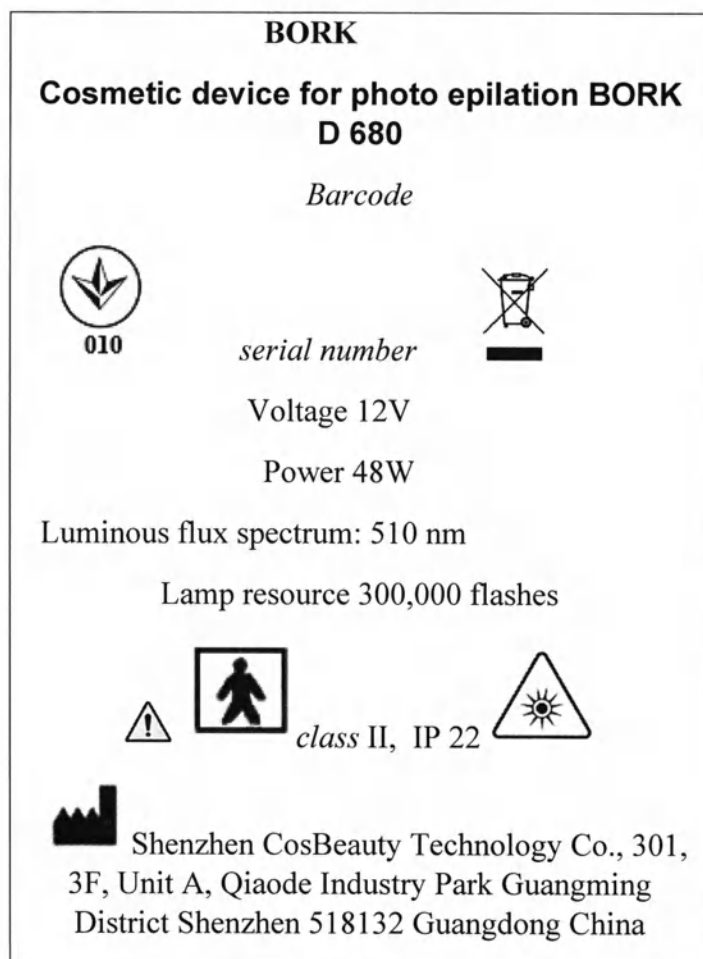1. Marking on the product packaging

BORK

INSTRUCTION FOR USE

Cosmetic device for photo epilation BORK D680

Cosmetic device for photo epilation BORK D 680

Barcode

Serial number

Adapter

Input voltage: 100-240V, output voltage: 12V

Apparatus

Voltage 12V

Power 48W

Luminous flux spectrum: 510 nm

Lamp resource 300,000 flashes

Shenzhen CosBeauty Technology Co., 301, 3F, Unit A, Qiaode Industry Park
Guangming District Shenzhen 518132 Guangdong China

010

class II, IP 22

Registration certificate: No. P3H _____ dated _____

Authorized representative, importer: LLC "BORK-Import", 108811, Moscow, settlement
Moskovsky, block 32, possession 17a, building 1, floor 3, room 104

Tel.: 8-800-500-88-99.

Transportation and storage conditions:

Temperature 5 °C - 30 °C

Relative humidity Not more than 65%

It is recommended to store the device in a cool dry place, isolated from places where acids and alkalis are stored.

Protect the device from significant temperature fluctuations and direct sunlight.

Nameplate for the AC adapter

BORK

INSTRUCTION FOR USE

Cosmetic device for photo epilation BORK D680

AC Adapter

100-240W

50/60Hz

1,2 A

12V

Made in China

Model K48V120400G white

Symbols description

Symbol	Description
	Caution! Refer to instructions for use
	Type BF
	Optical radiation
	Manufacturer
	Caution! Special disposal required!
	National product conformity mark (Ukraine)
	Class II equipment
	Indoor installation equipment
	TUV SVD certified GS certification

INSTRUCTION FOR USE

Cosmetic device for photo epilation BORK D680

	Double Insulation
---	--------------------------

PACKAGE

The packaging of the system provides protection against mechanical factors during transportation and storage, as well as the fullest use of the carrying capacity (capacity) of vehicles and the convenience of loading and unloading operations.

The device is packed in a protective chlorinated polyethylene bag (CPE bag). Consumer packaging is a cardboard box in which the elements of the device are placed.

The shipping container is a box made of corrected cardboard, designed for transportation by any mobile transport.

Number of pieces in a transport package 12 pieces

Make sure there are no visible signs of damage to the device or the contents of the package. If some of the contents are missing or there are visible signs of damage, contact an authorized representative of the manufacturer.

EMS

Guidance and MANUFACTURER'S declaration – ELECTROMAGNETIC EMISSIONS –

for all ME EQUIPMENT and ME SYSTEMS

Guidance and manufacturer's declaration – electromagnetic emissions		
The Cosmetic device for photo epilation BORK D680 is intended for use in the electromagnetic environment specified below. The customer or the user of the Cosmetic device for photo epilation BORK D680 should assure that it is used in such an environment.		
Emissions test	Compliance	Electromagnetic environment – guidance
RF emissions CISPR 11	Group 1	The Cosmetic device for photo epilation BORK D680 uses RF energy only for its internal function. Therefore, its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment

INSTRUCTION FOR USE

Cosmetic device for photo epilation BORK D680

RF emissions CISPR 11	Class B	The Cosmetic device for photo epilation BORK D680 is suitable for use in all establishments including domestic establishments and those directly connected to the public low-voltage power supply network that supplies buildings used for domestic purposes.
Harmonic emissions IEC	Class A	
Voltage fluctuations/flick er emissions IEC 61000-3-3	Complies	

INSTRUCTION FOR USE**Cosmetic device for photo epilation BORK D680****Guidance and MANUFACTURER'S declaration – electromagnetic IMMUNITY –
for all ME EQUIPMENT and ME SYSTEMS****Guidance and manufacturer's declaration – electromagnetic immunity**

The Cosmetic device for photo epilation BORK D680 is intended for use in the electromagnetic environment specified below. The customer or the user of The Cosmetic device for photo epilation BORK D680 should assure that it is used in such an environment.

Immunity test	IEC 60601 test level	Compliance level	Electromagnetic environment – guidance
Electrostatic discharge (ESD) IEC 61000-4-2	± 6 kV contact ± 8 kV air	± 6 kV contact ± 8 kV air	Floors should be wood, concrete or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30 %.
Electrical fast transient/burst IEC 61000-4-4	± 2 kV for power supply lines ±1 kV for input/output lines	± 2 kV for power supply lines ±1 kV for input/output lines	Mains power quality should be that of a typical commercial or hospital environment.
Surge IEC 61000-4-5	±1 kV different mode ±2 kV common mode	±1 kV different mode ±2 kV common mode	Mains power quality should be that of a typical commercial or hospital environment.
Voltage dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11	<5% UT; (>95% dip in UT) for 0.5cycle <40% UT; (>60% dip in UT) for 5cycles <70% UT; (>30% dip in UT) for 25cycles <5% UT; (>95% dip in UT) for 5	<5% UT; (>95% dip in UT) for 0.5cycle <40% UT; (>60% dip in UT) for 5cycles <70% UT; (>30% dip in UT) for 25cycles <5% UT; (>95% dip in UT) for 5	Mains power quality should be that of a typical commercial or hospital environment. If the user of The Cosmetic device for photo epilation BORK D680 requires continued operation during power mains interruptions, it is recommended that the Cosmetic device for photo epilation BORK D680 be

INSTRUCTION FOR USE

Cosmetic device for photo epilation BORK D680

	sec	sec	powered from an uninterruptible power supply or a battery
Power frequency (50/60 Hz) magnetic field IEC 61000-4-8	3 A/m	3 A/m	Power frequency magnetic fields should be at levels characteristic of a typical location in a typical commercial or hospital environment.

NOTE - UT is the AC mains voltage prior to application of the test level.

INSTRUCTION FOR USE

Cosmetic device for photo epilation BORK D680

Guidance and MANUFACTURER'S declaration – electromagnetic IMMUNITY –
for ME EQUIPMENT and ME SYSTEMS that are not LIFE-SUPPORTING

Guidance and manufacturer's declaration – electromagnetic immunity

The Cosmetic device for photo epilation BORK D680 is intended for use in the electromagnetic environment specified below. The customer or the user of the V30 should assure that it is used in such an environment.

IMMUNITY test	IEC 60601 TEST LEVEL	Compliance level	Electromagnetic environment – guidance
Conducted RF IEC 61000-4-6	3V rms 150 kHz to 80 MHz	3V	Portable and mobile RF communications equipment should be used no closer to any part of the V30, including cables, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter. Recommended separation distance $d = \left[\frac{3.5}{V_1} \right] \sqrt{P}$
Radiated RF IEC 61000-4-3	3V/m 80 MHz to 2.5GHz	3V/m	$d = \left[\frac{3.5}{E_1} \right] \sqrt{P} \quad 80 \text{ MHz to } 800 \text{ MHz}$ $d = \left[\frac{7}{E_1} \right] \sqrt{P} \quad 800 \text{ MHz to } 2.5 \text{ GHz}$ where p is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer and d is the recommended separation distance in metres (m)^b. Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey^a, should be less than the compliance level in each frequency range^b. Interference may occur in the vicinity of equipment marked with the following symbol: 

NOTE 1 At 80 MHz and 800 MHz, the higher frequency range applies.

INSTRUCTION FOR USE

Cosmetic device for photo epilation BORK D680

NOTE 2 These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people

a Field strengths from fixed transmitters, such as base stations for radio (cellular/cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast and TV broadcast cannot be predicted theoretically with accuracy. To access the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the location in which the Cosmetic device for photo epilation BORK D680 is used exceeds the applicable RF compliance level above, the Cosmetic device for photo epilation BORK D680 should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as reorienting or relocating the Cosmetic device for photo epilation BORK D680.

b Over the frequency range 150kHz to 80MHz, field strengths should be less than 3V/m

Recommended separation distances between portable and mobile RF communications equipment and the ME EQUIPMENT or ME SYSTEM – for ME EQUIPMENT and ME SYSTEMS that are not LIFE-SUPPORTING

Recommended separation distances between

portable and mobile RF communications equipment and the The Cosmetic device for photo

The Cosmetic device for photo epilation BORK D680 is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. The customer or the user of the The Cosmetic device for photo epilation BORK D680 can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters) and The Cosmetic device for photo epilation BORK D680 as recommended below, according to the maximum output power of the communications equipment

Rated maximum output power of transmitter W	Separation distance according to frequency of transmitter m		
	150 kHz to 80 MHz $d = [\frac{3.5}{V_1}] \sqrt{P}$	80 MHz to 800 MHz $d = [\frac{3.5}{E_1}] \sqrt{P}$	800 MHz to 2,5 GHz $d = [\frac{7}{E_1}] \sqrt{P}$
0.01	0.12	0.12	0.23
0.1	0.38	0.38	0.73
1	1.2	1.2	2.3
10	3.8	3.8	7.3
100	12	12	23

For transmitters rated at a maximum output power not listed above, the recommended separation distance d in metres (m) can be estimated using the equation applicable to the frequency of the transmitter, where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer.

NOTE 1 At 80 MHz and 800 MHz, the separation distance for the higher frequency range applies.

NOTE 2 These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people

PROCEDURE AND CONDITIONS FOR DESTRUCTION OR UTILIZATION OF MEDICAL DEVICE

If a product reaches the end of its useful life, the user must send the product to distributors or dispose of according to the national electronics label for disposal.

NAME AND LEGAL ADDRESS OF MANUFACTURER AND AUTHORIZED

REPRESENTATIVE IN RUSSIAN FEDERATION

MANUFACTURER:

Shenzhen CosBeauty Technology Co., Ltd., China

301, 3F, Unit A, Qiaode Industry Park Guangming District Shenzhen 518132 Guangdong China

Manufacturing site:

Shenzhen CosBeauty Technology Co., Ltd., China

301, 3F, Unit A, Qiaode Industry Park Guangming District Shenzhen 518132 Guangdong China

AUTHORIZED REPRESENTATIVE IN RUSSIAN FEDERATION:

Limited Liability Company "BORK-Import"

108811, Moscow city, Moskovsky settlement, block 32, building 17a building 1, floor 3 room 104

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