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CERTIFICATE

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China Council for the Promotion of International Trade
China Chamber of International Commerce

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THIS IS TO CERTIFY THAT: the seal of SHENZHEN FANSIZHE SCIENCE AND TECHNOLOGY CO., LTD on the annexed DOCUMENT is genuine.

China Council for the Promotion
of International Trade

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Authorized
Signature:

Gao Lichao

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(Date: Jul. 19, 2024)

证 询网址 Website for verifying the certificate: <http://www.rzccpit.com/validate.html>

Beauty requires
no sacrifice!

If you experience difficulty using our devices, please call the Kitfort Hotline before contacting the shop:

8-800-775-56-87

(Mon-Fri from 9:30 a.m. to 5:30 p.m. Moscow time).

info@kitfort.ru

We will explain the peculiarities of the appliance's operation and advise you on any other issues

Contents

Information	4
Set	4
Ind area of application	4
Uses for use	5
Indications to use	5
Possible complications	5
Copilot design	8
Preparation and use	10
Cleaning and maintenance	14
Care and storage	14
Troubleshooting	14
Technical data	16
Precautions	16
Manufacturer's information	18
Manufacturer's warranty	18

General information

Jeffort KT-3129 photoepilator (hereinafter photoepilator) makes it painless and effective to get rid of unwanted hair for a long time, and you do not have to spend time visiting a beauty salon.

For professional hair removal at home, the photoepilator uses Intense Pulsed Light (IPL) technology. The operating principle is based on the effect of light pulse on melanin in the hair. Melanin absorbs the light pulse and converts it into heat. The heat goes to the root of the hair, where it heats up and destroys the hair follicle. When using IPL technology, the result is achieved only when removing dark hair on the skin from light to medium skin tones. Photoepilators using this technology are not suitable for the removal of blonde, red and grey hair, or for use on very dark skin.

Photoepilation will not only get rid of unwanted hair, but will also improve skin elasticity and blood microcirculation, as well as reduce wrinkles by stimulating collagen production during the procedure.

If you undergo a course of photoepilation, you can achieve smooth skin for a long time. Unlike epilators, shaving machines, sugaring and waxing, photoepilation does not damage the skin and, if done correctly, does not cause any skin irritation.

You can choose one of 8 levels of intensity of exposure, focusing on the recommendations in the manual and your own feelings. The device has two modes of operation: manually operated and automatic. The manually operated mode is designed to treat small areas of skin: bikini zone, underarms, knees, and the automatic mode is suitable for treating large areas, such as arms and legs.

The photoepilator is equipped with accidental switch-on protection: the device will not produce a light pulse if the sensor built into the photoepilator is not in contact with the skin. Special glasses protect your eyes from the light pulses during the procedure.

Complete set

1. Photoepilator - 1 pc.
2. Power supply unit - 1 pc.
3. User manual - 1 pc.

Accessories:

4. Safety glasses - 1 pc.
5. Shaving stick with cover - 1 pc.

Purpose and area of application

Photoepilator KT-3129 (hereinafter - photoepilator) is an electronic device designed to remove unwanted hair and intended for use without a doctor's prescription.
Area of application - for use at home.

Indications for use

- ingrown hair problem;
- hypertrichosis (excessive hair on different parts of the body);
- male pattern hair distribution in women.

Contraindications for use

- presence of tattoos in the treatment area;
- freckles in the treatment area;
- birthmarks in the treatment area;
- moles and/or warts;
- thyroid disorders;
- pregnancy;
- lactation;
- fresh tan;
- abrasions and cuts;
- mechanical injuries;
- varicose veins;
- psoriasis;
- eczema;
- diabetes mellitus;
- hypertension;
- cardiovascular disease;
- wearing of a pacemaker;
- infectious and viral diseases;
- epilepsy;
- blood coagulation disorders;
- tendency to manifest keloid and hypertrophic scars;
- mental disorders;
- malign tumors;
- hypersensitivity of the skin to light;
- taking anti-inflammatory, diuretic, hypotensive drugs;
- treatment with medications that contain iron and magnesium;
- oncology, including within 3 months after chemotherapy.

Possible complications

- redness of the skin on the treated skin area;
- a burning sensation after the treatment;
- possible development of folliculitis (inflammation of hair follicles);
- temporary hypo- or hyperpigmentation of the treated skin areas;
- eye damage and reduced visual acuity, occur if procedures are carried out without protective goggles;
- ingrown hairs in the epilation area.
- For 3 days after photoepilation, it is not allowed to visit a bathhouse, sauna, swimming pool, and it is recommended to limit water procedures as much as possible.

ATTENTION!

DO NOT look into the light window on the body of the photoepilator when the power is on and during its use.

DO NOT direct the light emitted by the photoepilator into your eyes.

ONLY activate the flash when the photoepilator is pressed firmly against the skin area to be treated. Turn away when the photoepilator is pointed at the skin.

DO NOT treat the skin around the eyes (eyebrows or eyelashes).

The photoepilator can be used to treat the upper lip, chin, cheeks, jaw and neck in women.

DO NOT use the photoepilator to epilate surfaces other than skin.

Safety precautions

The device is intended for individual use.

It is not intended for use by persons under 18 years of age.

DO NOT USE the device if the skin in the desired treatment area has been peeled.

Treatment of tanned skin with Intense Pulsed Light (IPL) may result in permanent hyperpigmentation (darkening) or hypopigmentation (lightening) of the treated skin.

DO NOT USE the device if you know your skin is sensitive to sunlight (photodermatosis) or you are taking medication that makes your skin more sensitive and causes photosensitisation, such as, Retin A, Accutane and/or other retinoids for external use.

device on photosensitive skin may cause damage such as swelling or blistering.
The device is not intended for use by persons with reduced physical, sensory or mental capabilities unless supervised or have been given advice on safety and proper use of the device by a person responsible for safety.

Ensure that children do not use the device as a toy. Keep the device out of the reach of children at all times.
Check the device and the cords regularly for visible damage.

If the device is damaged or cracked, stop using it.

Using a damaged device may result in an electric shock. If there are any defects, take the device to an authorised service centre for inspection and repair.

Repair of the device by a non-specialist may result in serious risk to the user. Always unplug the device after use.

WARNING!

Never attempt to repair or disassemble the device yourself, as it contains dangerous electrical parts and pulsed light energy, contact with which may cause serious injury to the body and or severe eye damage.

Do not use in, on, or near wet surfaces, such as in a bathtub, shower, sink, or other receptacles that may contain water. Do not use a device into which water has penetrated. Electric shocks may occur if you use a device into which water has penetrated.

Do not use a damaged device (e.g., cracked body, damaged cord, etc.). Internal parts contain an electrical charge that may be dangerous.

Do not use if the glass light filter is cracked or missing.

Check the glass light filter regularly for damage.

Do not touch the glass filter immediately after using the device, as it will heat up during use.

Do not treat the device with cleaning agents: petroleum-based cleaning agents or flammable cleaning agents (alcohol, acetone) as they are flammable.


Do not use abrasive cleaning agents. They may damage the device or leave scratches on the glass light filter.

Do not attempt to use the device for purposes other than removing unwanted hair.

Do not use the device on surfaces other than skin. Do not press the treatment button if the treatment surface does not fully adhere to the skin.

Do not block the cooling holes. The device may overheat and pose a risk of electric shock or fire.

Photoepilator design

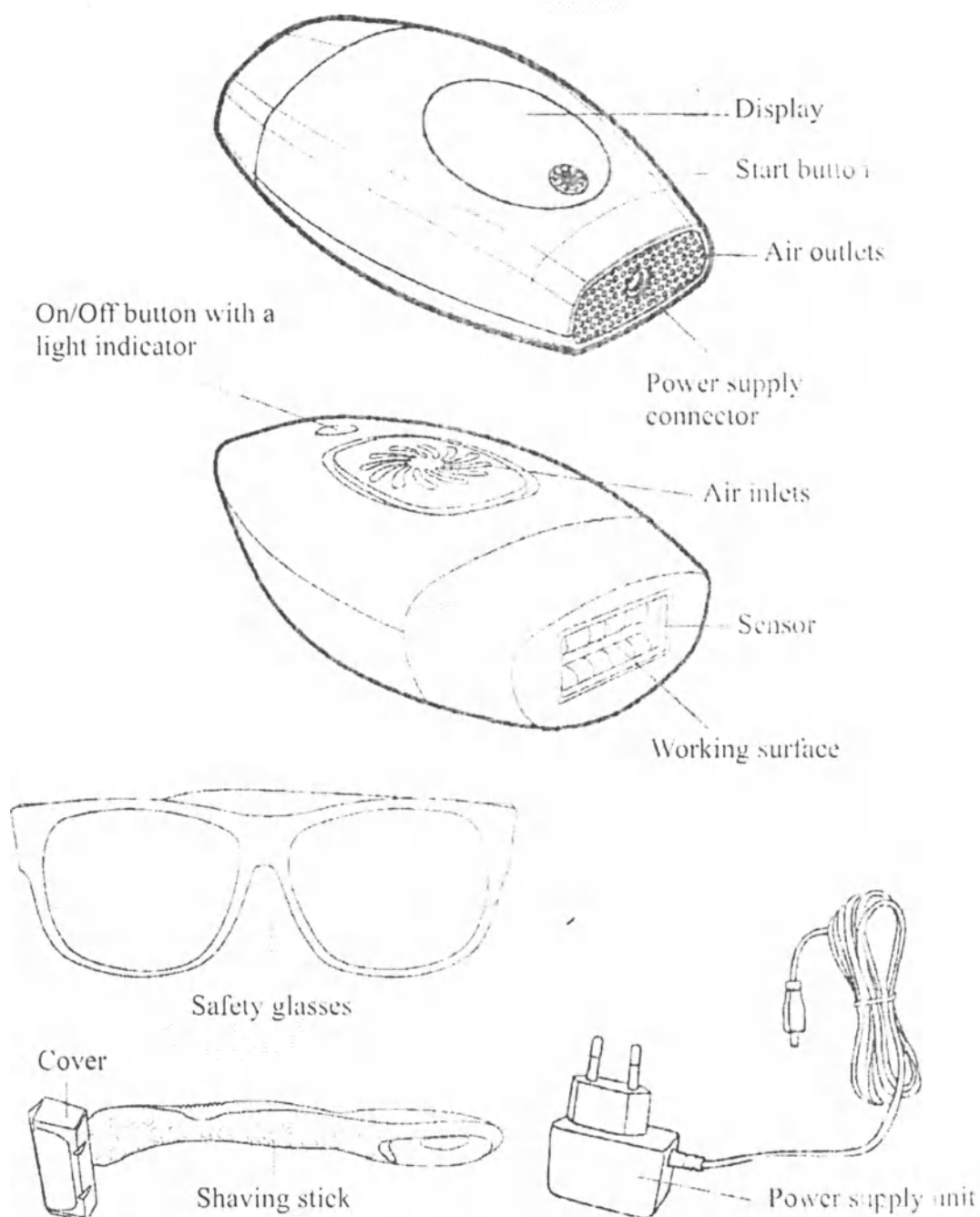The operating principle of a photoepilator is quite simple. When the device is switched on, the built-in broadband lamp sends a powerful light pulse (flash), which affects the melanin in the hair. The melanin absorbs the light pulse and converts it into heat. The heat travels to the root of the hair, where it heats up and destroys the hair follicle.

the hair, the more melanin it contains. Hair containing a high amount of melanin absorbs more the light pulse, so the hair follicle is destroyed more easily. The procedure is not effective in red, blonde and grey hair, as they contain little or no melanin. The skin in the treatment area is not affected because the skin contains little or no melanin (except for very dark skin).


During photoepilation, safety glasses will protect your eyes from the harmful effects of residual light from flashes.

Before photoepilation you can remove hair with the shaving stick from the set or use your own.

Air enters the unit through the air inlets, cools the internal parts of the unit and exits through the holes on the bottom of the unit.

To avoid accidents, the light pulses (flashes) only when the sensor is pressed against the skin area to be treated. The sensor is located along the perimeter of the working surface.

Display

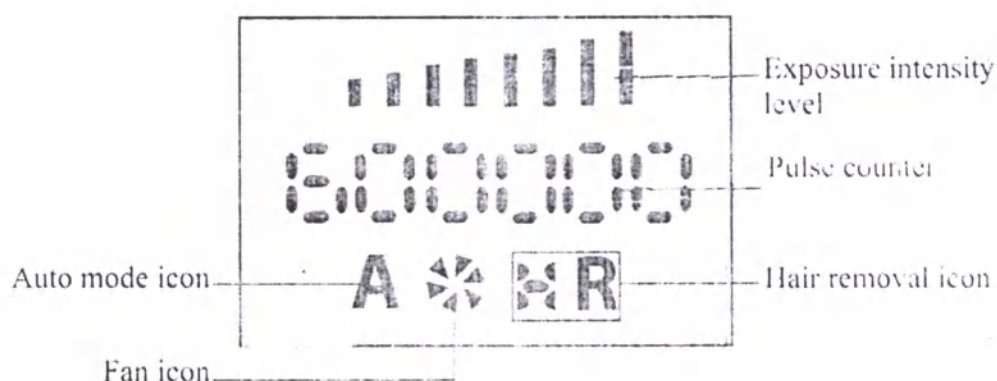The *display* shows the intensity level, pulse counter, and auto mode, hair removal and fan icons.

A long press on the **On/Off button** switches the photoepilator on and off and a quick single press switches the intensity level from 1 to 8. The intensity of the treatment is determined by the amount of radiation energy. The higher the intensity level, the greater the radiation energy and the stronger the effect. The intensity level corresponds to the number of divisions on the device display. The hair removal and fan icons light up on the display when the device is switched on.

Selecting the power level: there are 8 power levels according to different skin types, 1 being the lowest and 8 being the highest. As darker skin absorbs energy more actively, we recommend using the lowest levels for this skin type. We suggest setting the levels in accordance with Table 1.

ive increase of the intensity level may cause burns, so use the intensity level that suits your skin (see

Using the **Start button** activates the pulse in manually operated mode and cycles through the operating modes. The **Auto** icon appears when **Auto mode** is selected.

The **pulse counter** shows the number of remaining flashes the lamp can still produce. The lamp is designed to produce 600,000 flashes. For one course of epilation, consisting of hair removal from arms, legs, bikini area and underarms, about 600 flashes on average are required (the number of flashes is individual for each person and depends on the type and density of hair, as well as the treatment area). That is, the lamp's lifespan is enough to carry out 1000 procedures of photoepilation. Assuming that a full course of photoepilation requires 8 treatments, the lamp's lifespan will be enough for 125 courses.

Note. When the lamp's lifespan is exhausted and the error code "ERR0" appears on the display, stop using the device and dispose of it in the safest way possible, e.g. by taking it to a special recycling centre.

Operating modes

The **manual mode** is suitable for small areas of skin: bikini area, underarms, knees.

Once the working surface of the device with the sensor is pressed against the skin, the indicator light built into the On/Off button will start flashing - the device is ready to produce a pulse. If the indicator light is on continuously, the sensor is not pressed firmly enough against the skin.

To produce a single light pulse, press the Start button. After the pulse is delivered, the indicator light will go out and light up again after 3 seconds, indicating that the device is ready for use again.

Move the working surface of the photoepilator to the next skin area and activate the single light pulse again. In this manner, treat the entire selected area.

The automatic mode is convenient for removing hair from larger areas such as arms and legs.

When the sensor comes into contact with the skin, the photoepilator starts to flash automatically and the indicator light is not illuminated. Slowly move the working surface of the photoepilator over the skin from bottom to top. In this way, treat the entirety of the chosen area.

Preparation and use

Preparing for use

Unpack the photoepilator, power supply unit, safety glasses and shaving stick from the box and remove all packaging materials.

Shave the hair 1-2 days before the photoepilation. For effective photoepilation, hair should be 1-2 mm long.

is not visible, there will be no effect from the procedure, as the light energy will not be able to reach the hair. If the hair is longer than necessary, the energy of the light pulse will dissipate before it reaches the hair and will not be able to destroy it.

Application precautions

If you are not sure if you are allowed to have photoepilation, consult a dermatologist.

To avoid the harmful effects of residual light from the flashes on your eyes, only use the device wearing protective glasses and do not use the photoepilator with other people or pets nearby.

One month before photoepilation and during the course of photoepilation it is forbidden to remove hair with epilator, have sugaring procedures or wax the skin.

A month before photoepilation and for 2 weeks after, do not expose your skin to direct sunlight, do not visit a solarium, do not peel or scrub your skin.

During the course of photoepilation exclude from your diet foods that may cause photoallergic reactions (e.g., celery, parsley, fennel, figs, some citrus fruits), do not allow contact with plants whose juice can increase photosensitivity (eg, hogweed).

Do not perform photoepilation on minors.

Do not use the photoepilator on Fitzpatrick phototype VI skin.

Do not use the device:

- on areas with tattoos, freckles, birthmarks, moles or warts;
- on areas with eczema, psoriasis, abrasions, open wounds or infections (herpes), abnormal skin conditions caused by systemic or metabolic diseases (e.g. diabetes);
- during pregnancy and lactation;
- if you have epilepsy, impaired collagen production or blood clotting, immunodeficiency;
- if you have diseases associated with sensitivity to light and photosensitivity (e.g. porphyria disease, polymorphic photodermatitis, solar urticaria, lupus, etc.);
- if you have oncology, including within 3 months of receiving chemotherapy;
- if you have a pacemaker, neurostimulator, etc.;
- if you are taking medications that affect sensitivity to light and photosensitivity (e.g. antibiotics and hormones).

For 3 days after photoepilation, it is not allowed to visit a bathhouse, sauna, swimming pool, and it is recommended to limit water procedures as much as possible.

Before use

Before the treatment, carry out a test to check your skin's sensitivity to photoepilation. To do this, treat a small area of skin in manually operated mode at the first intensity level as described in the chapter "Use". If you do not feel any discomfort, treat other areas of your skin, gradually increasing the intensity level until you reach the highest level that is comfortable for you. After 48 hours, examine the skin where the test was performed. If the skin looks normal, you can use the photoepilator. However, you can only re-treat the area on which the test was carried out one week after the test.

Use

1. Make sure that the skin on which you are going to perform photoepilation is clean and dry. There should not be any cosmetic products on the skin.
2. Clean the working surface of the device with a medical alcohol wipe or a cotton disc soaked in chlorhexidine.
3. Connect the power unit to the power supply connector on the photoepilator and plug the power unit into the mains. The indicator light will start flashing.
4. Put on safety glasses.
5. Press and hold the On/Off button to switch the unit on. The display will illuminate. The indicator light will be on continuously. You will hear the sound of the fan running. By default, the first intensity level and manually operated mode are selected.
6. Press the On/Off button to change the intensity level. When selecting the intensity level, use the recommendations in the table in the Advice chapter.
7. Press and hold the Start button until the Auto mode icon lights up if you want to go into Auto mode. To return to manually operated mode, press the Start button once more. The Auto mode icon will disappear.
8. Press the working surface of the photoepilator against your skin at a 90° angle and treat the selected area of skin as described in the Operating modes chapter.

Warning! The device will not emit a light pulse if the sensor is not pressed against the skin, but the residual light from the flash can be harmful to the eyes. To prevent the eyes from being affected by the residual light from the flash, always wear safety glasses when using the device and do not use the photoepilator with other people or pets nearby.

Note. You may feel warmth and tingling during photoepilation and for a few minutes afterwards - that is normal. However, if the sensation is too painful, reduce the intensity or stop using the device and consult your doctor.

and hold the On/Off button to switch off the device. The display will turn off and the indicator
link.
plug the power supply unit from the mains. The indicator light will go out.

ce

average course of photoepilation is 6-10 treatments. Carry out the photoepilation procedure once every
ks for the first 6 times, then once every 4 weeks until you achieve the desired result. If necessary, after
pleting the course of photoepilation, maintain the result with 1 treatment every 1-2 months.

The duration of the course of photoepilation depends on the colour and thickness of hair, the amount of
in the growth stage at the time of photoepilation, as well as biological factors affecting hair growth.

Do not perform photoepilation on the same area of skin more than once a week. Frequent use will not result
faster results, but may increase the risk of skin damage.

After the treatment, apply a regenerating cream (e.g. panthenol) to the treated skin area to avoid dry and
ky skin.

For 3 days after photoepilation it is forbidden to visit a bathhouse, sauna, swimming pool, and it is
commended to limit water procedures as much as possible.

Do not expose your skin to direct sunlight, do not visit a solarium, do not peel or scrub the skin 2 weeks
ter the procedure of photoepilation. Please note that at the beginning of the course of photoepilation a
cedure is carried out every fortnight, so for the first 3 months you should not expose the skin to direct
nlight, visit a solarium, peel or scrub the skin.

During the course of photoepilation it is prohibited to remove hair with epilator, have sugaring procedures
wax the skin.

The first 2 days after the procedure do not apply cosmetics containing alcohol or greasy creams to the
ated skin.

Do not use decorative cosmetics for 2-3 days after photoepilation on the face.

Slight redness of the skin after the procedure is normal. The redness subsides within 48 hours. Consult your
ector if the redness does not go away after 2-3 days.

You may feel warmth and tingling during photoepilation and for a few minutes afterwards – that is normal.

Table 1 gives recommendations for the procedure for different skin types. Use on any area of skin that is
rker in colour or relatively more pigmented may cause unwanted phenomena such as burns, blisters and
in discolouration (hyper- or hypopigmentation).

area	Skin type on Fitzpatrick scale	I and II	III	IV	V
	Exposure intensity level	6-8	5-7	4-6	1-4
	Number of pulses	3-4	2-3	2-3	2-3
legs	Exposure intensity level	6-8	5-7	1-4	1-4
	Number of pulses	3-4	2-3	2-3	2-3
e and neck	Exposure intensity level	6-8	4-6	1-4	1-4
	Number of pulses	3-4	2-3	2-3	2-3
kini zone	Exposure intensity level	6-8	4-6	1-4	1-4
	Number of pulses	3-4	2-3	2-3	2-3

Cleaning and maintenance

Clean the photoepilator with a soft, damp cloth and then a dry cloth. Do not immerse the device in water or other liquids. Do not allow water to enter the body, controls, power supply connector or power supply unit. Wipe the work surface first with a damp cloth and then with a dry soft cloth. Clean dust from the air outlet and air inlet with a soft-tipped brush.

Care and storage

Store the device out of direct sunlight in a dry and cool place out of the reach of children. When storing the photoepilator, safety glasses and power supply unit for long periods of time, place them in a box or bag to protect them from dust.

Troubleshooting

Device does not switch on, the indicator light is not on

Probable cause	Solution
Mains voltage	Check presence of voltage in mains
Power supply unit malfunction	Refer to a service centre

No light pulse is emitted

Probable cause	Solution
The sensor is not in contact with the skin	To avoid accidents, the light pulse is only emitted when the sensor is pressed against the skin area to be treated
The manual mode is selected	To produce a single light pulse, press the Start button. After the pulse, the indicator light will go out and come back on after 3 seconds, indicating that the unit is ready for use again.

Severe pain during use

Probable cause	Solution
The intensity level is set too high	Stop using the device. Wait until the skin has recovered (redness and irritation have subsided). The next time you use the device, select an intensity level lower than before

Insufficient effectiveness of the procedure

Probable cause	Solution
The intensity level is set too low	Increase the intensity level. But do not set the level too high to avoid damaging the skin
Hair too light, grey or red	The best results are achieved when removing dark hair
Hair length too short/long	For effective photoepilation, hair should be 1-2 mm long

and of fault occurs, the device will switch off and display an error message on the display or start. The following table will help you decipher the error code.

Error code	Probable cause	Solution
ERR0	Lamp life expired	The device is no longer suitable for use. Stop using the device and dispose of it in the safest way possible, e.g. by taking it to a special recycling centre.
ERR1	High voltage part error	Refer to a service centre
ERR2	Fan defective	Refer to a service centre
ERR3	Device overheated and switched off	Switch off the device and unplug it from the mains. Wait until the device has cooled down completely. If the error code persists, refer to a service centre
ERR4	Problem with the circuit board	Refer to a service centre

If your situation is not shown above, please email us at info@kitfort.ru with photos or videos of your problem. Please also send us a photo of the serial number sticker located on the bottom or back of the device body.

For questions about purchasing consumables or accessories, please contact us at osh@kitfort.ru.

Technical data

- Power supply unit:
 - Input: ~100-240 V, 50/60 Hz
 - Output: ~12 V, 3 A
 - Power supply: 36W
 - Protection class against electric shock:
- Power supply unit: II
- Photoepilator: III
- Number of operating modes: 2
- Number of intensity levels: 8
- Radiation energy: 1.5-3.9 J/cm²
- Area of epilation zone: 4.2 cm²
- Number of pulses: 600000
- Wavelength: 470-1100 nm.
- Noise level: 28-32 dB
- Cord length: 1.1 m
- Dimensions of the device: 84 x 50 x 143 mm
- Dimensions of packaging: 133 x 90 x 165 mm
- Net weight: 0.38kg
- Gross weight: 0.51kg

Precautions

Please read the operating instructions carefully. Pay particular attention to the safety precautions. Always keep the instruction manual handy.

- 1. The device is intended for use at home and can be used in flats, country houses, hotel rooms, offices and other similar places for non-industrial and non-commercial use.
- 2. Use the photoepilator only for its intended purpose and in accordance with the instructions in this manual. Any other use of the device will be considered a violation of proper operating practices.
- 3. Before connecting the device to an electrical outlet, make sure that the power supply parameters indicated on the device are the same as those of the power source being used.
- 4. To prevent electric shock, do not immerse the device in water or other liquids.
- 5. Do not carry the photoepilator by the power cord. Do not pull on the power cord when disconnecting the power supply unit from the device.
- 6. Do not use the device if the power cord, plug or other parts of the device are damaged. To avoid electric shock, do not disassemble the device by yourself – have it repaired by a qualified technician. Remember, incorrect assembly of the device increases the risk of electric shock during operation.

Persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge must use the photoepilator under the supervision of persons responsible for their safety or after being trained in the use of the device. Do not allow children to play with the device.

To avoid the harmful effects of residual light from the flashes on the eyes, do not use the photoepilator near people or pets nearby.

Do not leave the device unattended when it is in operation. Switch it off and unplug it from the mains power supply when not using the device for an extended period of time or before servicing. Remember, an unattended device can cause ignition!

1. Only place the photoepilator on a stable horizontal surface, not less than 10 cm away from the wall and edge of the table. To avoid blocking the air inlets and damaging the device, do not place the photoepilator on soft surfaces when turned on.

2. Do not attempt to bypass the device's power lock.

3. Do not drop the photoepilator or subject it to shocks.

4. Wear safety glasses when performing the procedure.

5. Keep the device out of the reach of children.

6. The device is not intended for use by children under 18 years of age.

7. One month before photoepilation and for 2 weeks after, do not expose your skin to direct sunlight, do not visit a solarium, do not peel or scrub the skin.

8. Do not perform photoepilation on skin phototype VI on the Fitzpatrick scale.

9. Do not use the device:

on areas with tattoos, freckles, birthmarks, moles or warts;

on areas with eczema, psoriasis, abrasions, open wounds or infections (herpes), abnormal skin conditions caused by systemic or metabolic diseases (e.g. diabetes);

during pregnancy and lactation;

if you have epilepsy, impaired collagen production or blood clotting, immunodeficiency;

if you have diseases associated with sensitivity to light and photosensitivity (e.g. porphyria disease, polymorphic photodermatitis, solar urticaria, lupus, etc.);

if you have oncology, including within 3 months after chemotherapy;

if you have a pacemaker, neurostimulator, etc.;

if you are taking medications that affect sensitivity to light and photosensitivity (e.g. antibiotics and hormones).

3 days after photoepilation it is forbidden to visit a bathhouse, sauna, swimming pool, and it is recommended to limit water procedures as much as possible.

Use only the accessories or components offered by the manufacturer. Use of other accessories may cause damage to the device or injury.

If the power supply unit is damaged, replace it with a special power supply unit or set obtained from the manufacturer or service department.

Service life: 2 years

Warranty period: 1 year

Manufacturer:

Shenzhen Fansizhe Science And Technology Co., Ltd.

Place of manufacture address: Second floor, third floor, fourth floor, No. 11, Henggangxia Street, Xinqiao Henggangxia Industrial Zone, Xiangshan Community, Xinqiao Street, Baoan District, Shenzhen

Authorised representative: Shenzhen AIFAN Technology Co., Ltd., 4 Floor, Huanan Industrial Park, Huananfeng 4 Road, Bashixia East, Fuyun Town, Baoan District, Shenzhen, China.

Authorised representative: AERO-TRADE LLC

OGRN 1099847029370, INN 7811451960, KPP 781301001

Address: 197022, Saint Petersburg, federal city intracity territory Aptekarsky Ostrov municipal district, Instrumentalnaya ul., 3, letter П, room 3-H

Hotline of the manufacturer's authorised representative:

8-800-775-56-87 (Mon-Fri from 9:30 to 17:30 Moscow time), info@kitfort.ru

You can get the addresses of service centres from the hotline operator or on the kitfort.ru website.

The month and year of manufacture are indicated on the underside of the packing box.

Warranty conditions

Mechanical damage to the body, accessories or parts of the device is not a warranty case.

Failure of the device as a result of water penetration into the device body, into the area of controls, into the area of the power supply connector or into the power supply unit is not a warranty case.

Lamp's lifespan is not a warranty case.

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