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

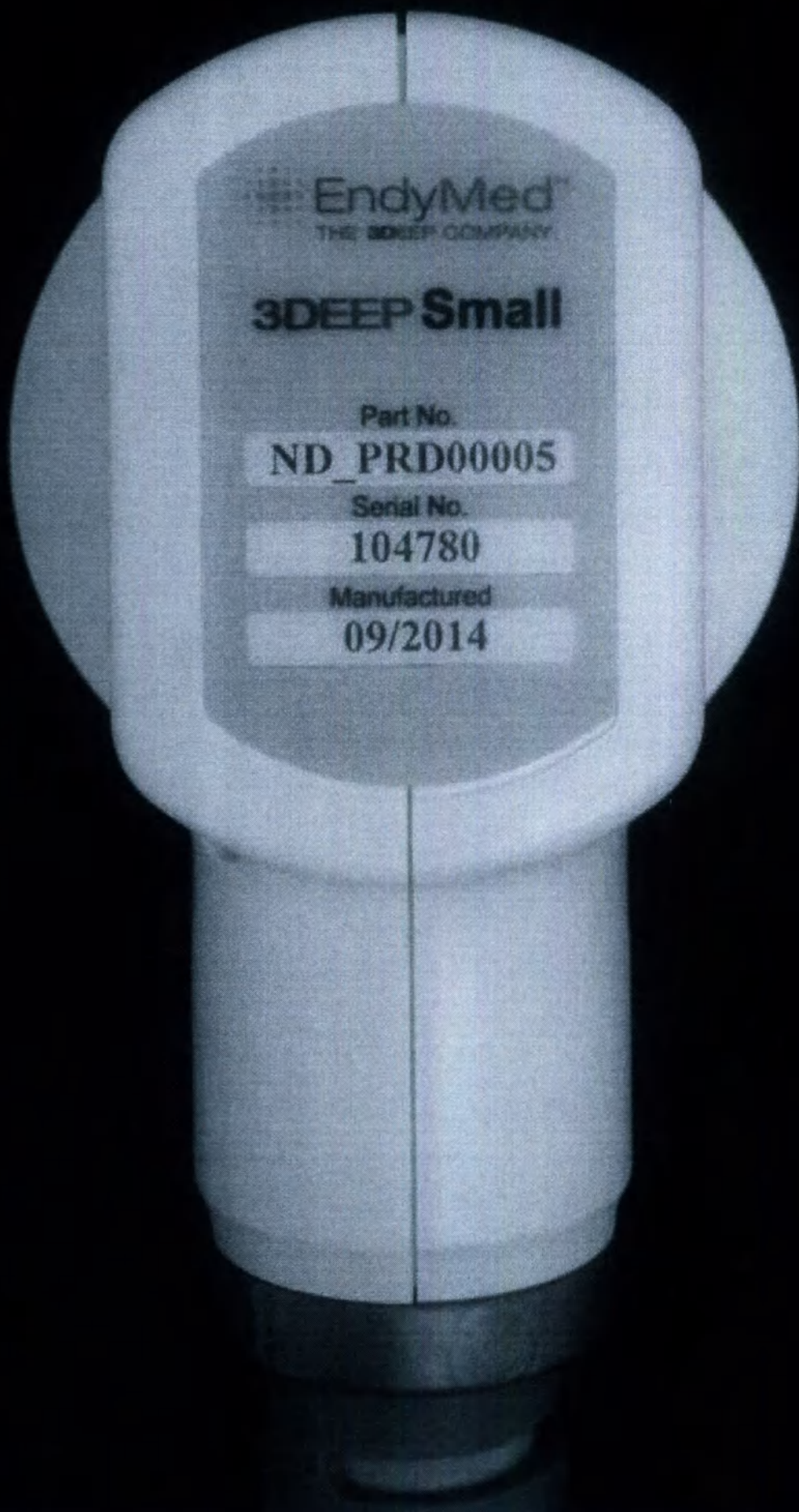
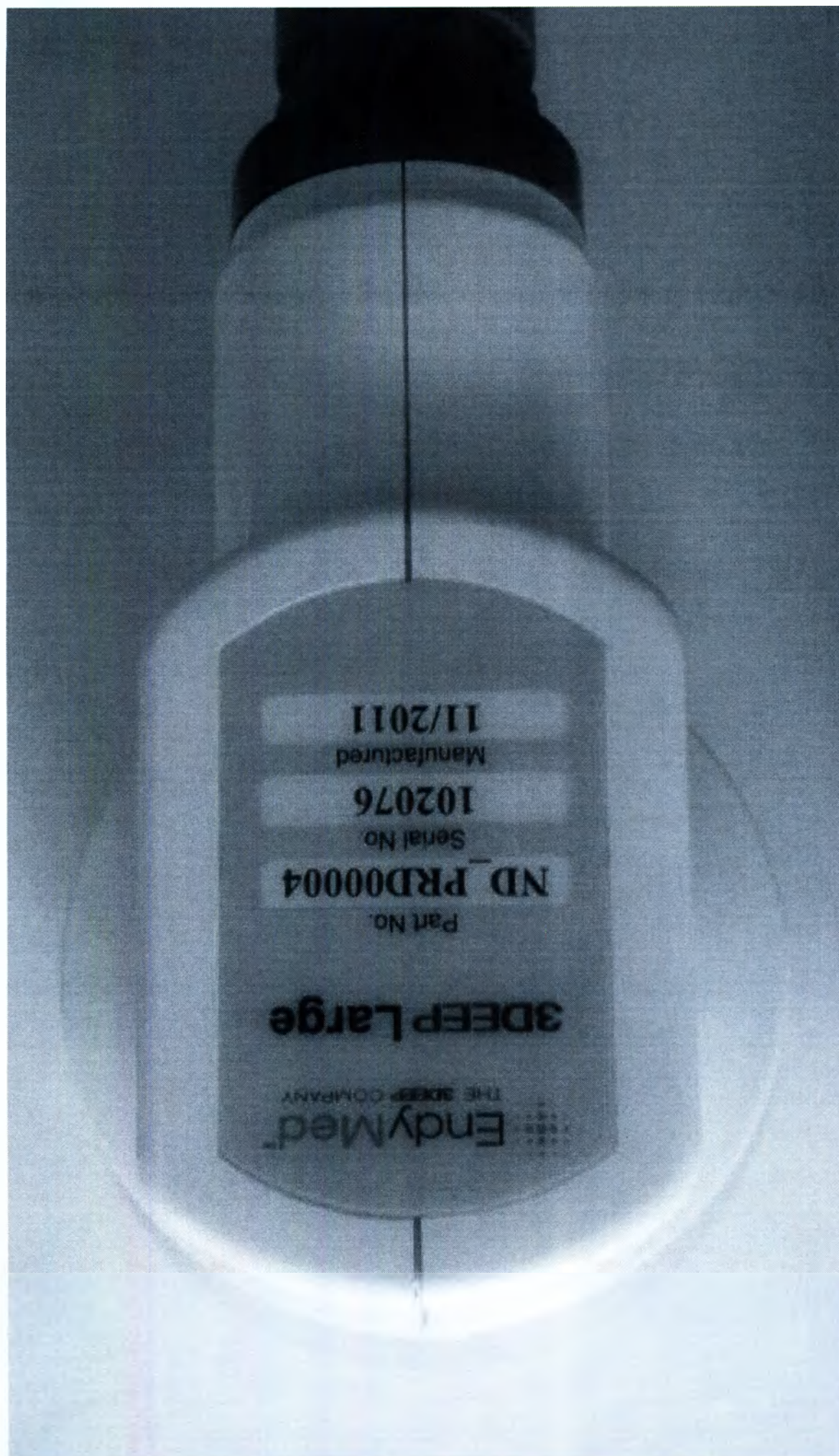
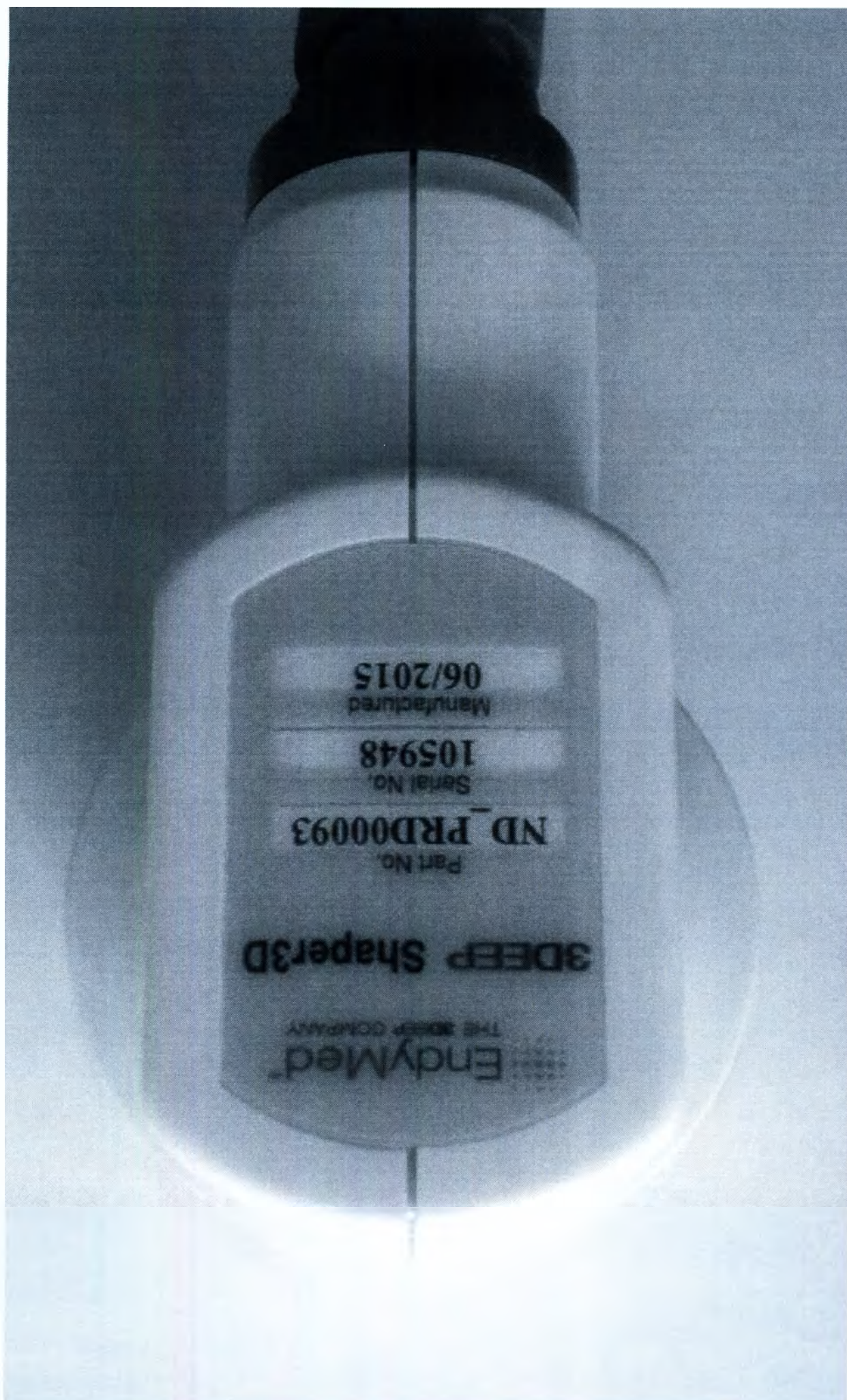
Manufactured

09/2016

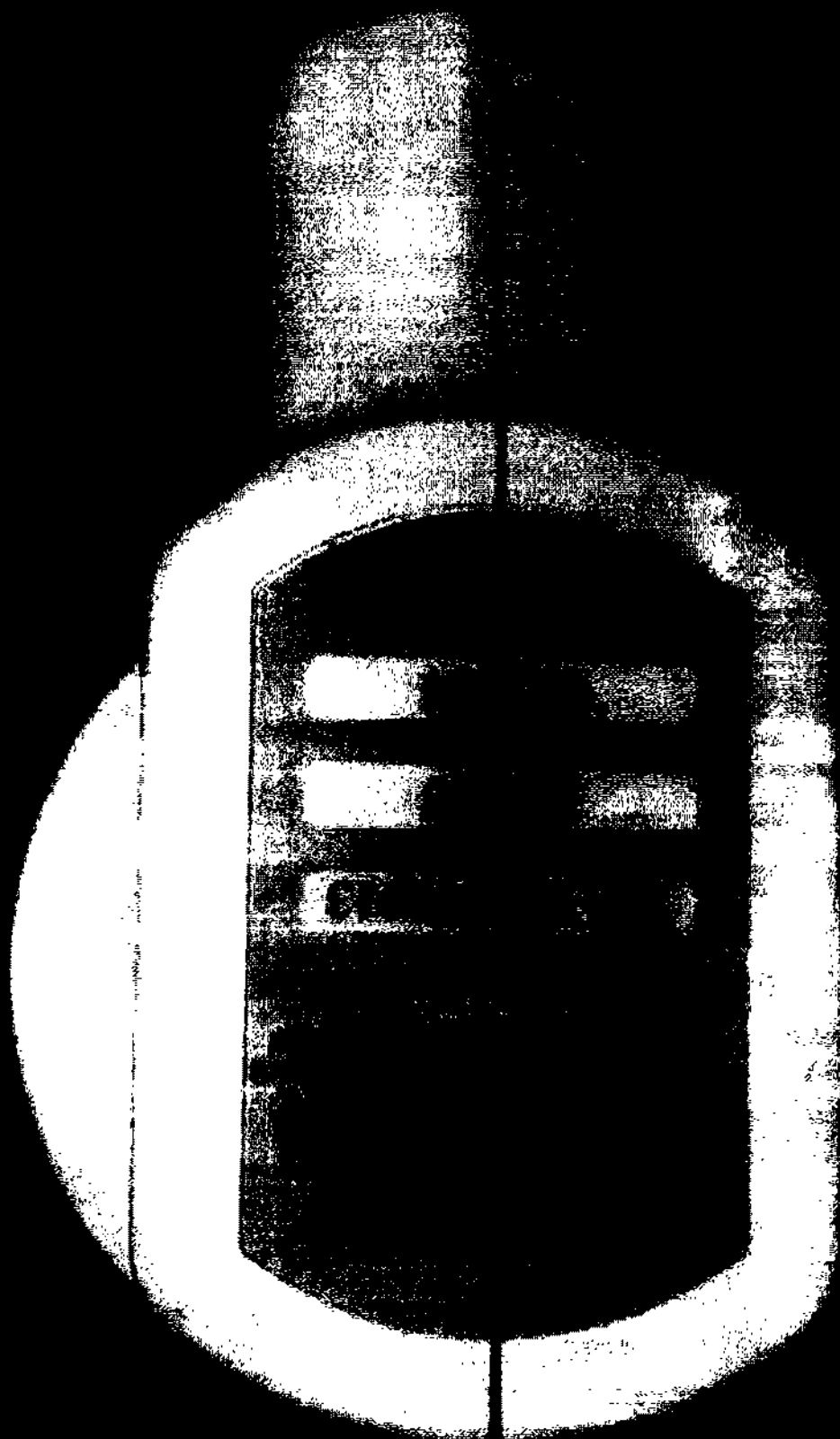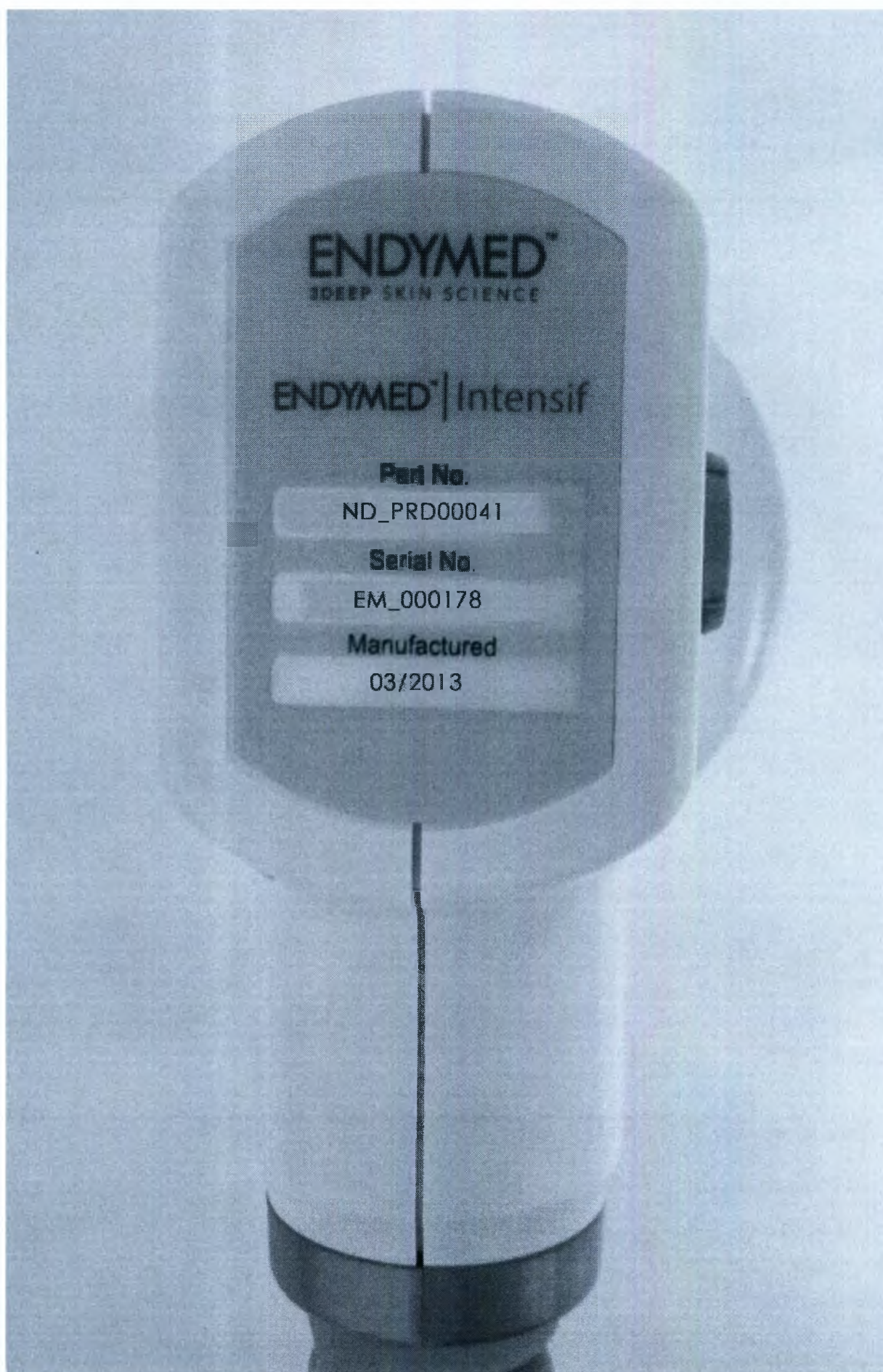10 Treatment Tips for the FSR Handpiece

up to

450

pulses per tip

EndyMed™
THE 3DEEP® COMPANY

10 FSR Treatment Templates

The FSR Template enables easy positioning of the FSR Handpiece during treatment to ensure uniform coverage of treatment area.

1. Place the Template on the standard treatment area.
2. Use a light, oil-free make-up base and a small sponge to spread the make-up through the Template holes, creating a uniform matrix grid on the skin.
3. Position one corner of the lower edge of the FSR tip on one of the outside dots of the grid, so the surface of the tip covers the first area to be treated. Deliver pulse.
4. Lift the FSR Handpiece and position the water tip corner on the next grid dot so that the lower edge of the tip covers the next area. Deliver pulse.
5. Repeat Step 4 for full and uniform coverage of treatment area.

Treatment Tip

Fractional Skin Resurfacing Handpiece

ENDYMED™

Please read the instructions carefully before use.

CE
0483

STERILE EO

12 Leshem St. Caesarea,
3066900, Israel

Do not use if
package is damaged.

02/2016

052016

02/2018

Treatment Tip

Fractional Skin Resurfacing Handpiece

ENDYMED™

Please read the instructions carefully before use.

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

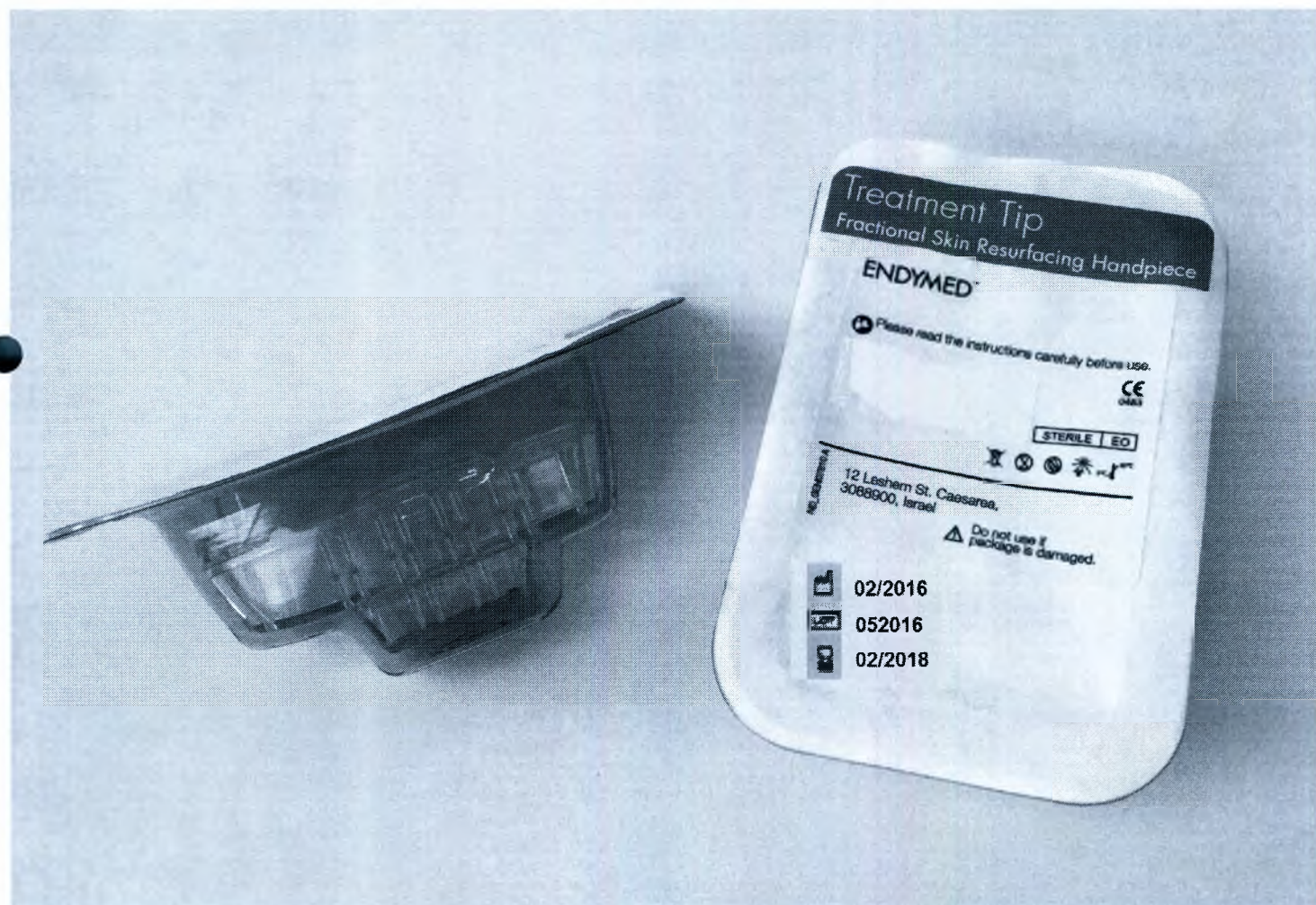
1. Place the Template on the intended treatment area.
2. Use a light, oil-free makeup base and a small sponge to spread the makeup through the template holes, creating a uniform matrix grid on the skin.
3. Position one corner of the lower edge of the FSR tip on one of the outside dots of the grid, so the surface of the tip covers the first area to be treated. Deliver pulse.
4. Lift the FSR Handpiece and position the same tip corner on the next grid dot so that the lower edge of the tip now covers the next area. Deliver pulse.
5. Repeat Steps 4 for all and obtain coverage of treatment area.

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
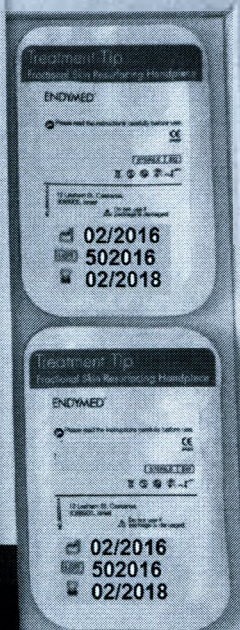
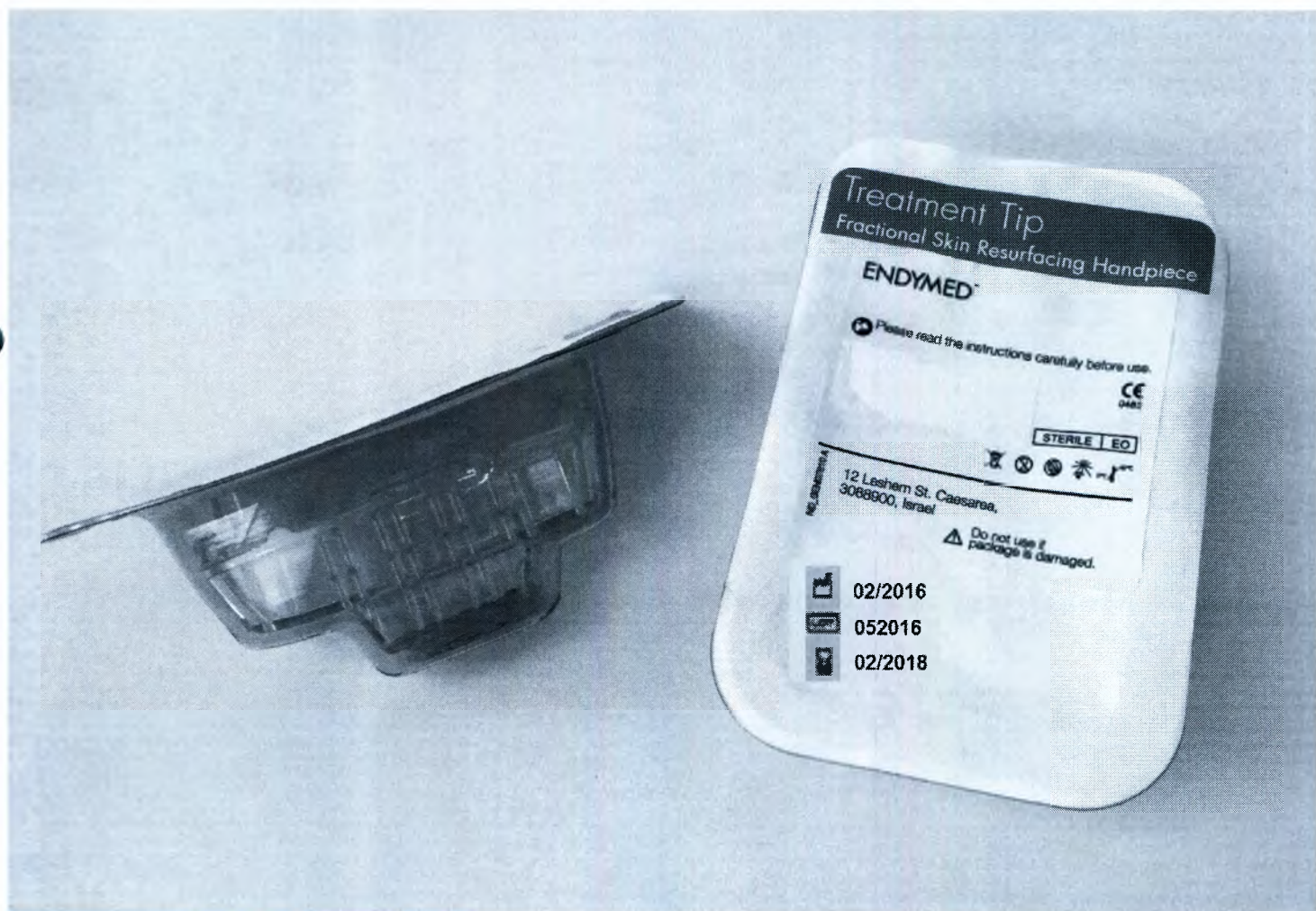
pulses per tip

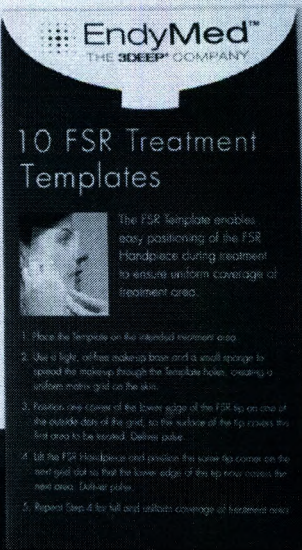
10 FSR Treatment Templates

The FSR template enables easy positioning of the FSR Handpiece during treatment to ensure uniform coverage of treatment area.

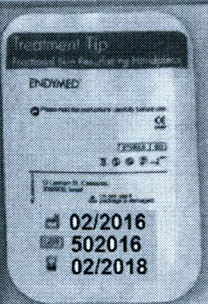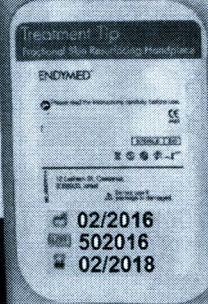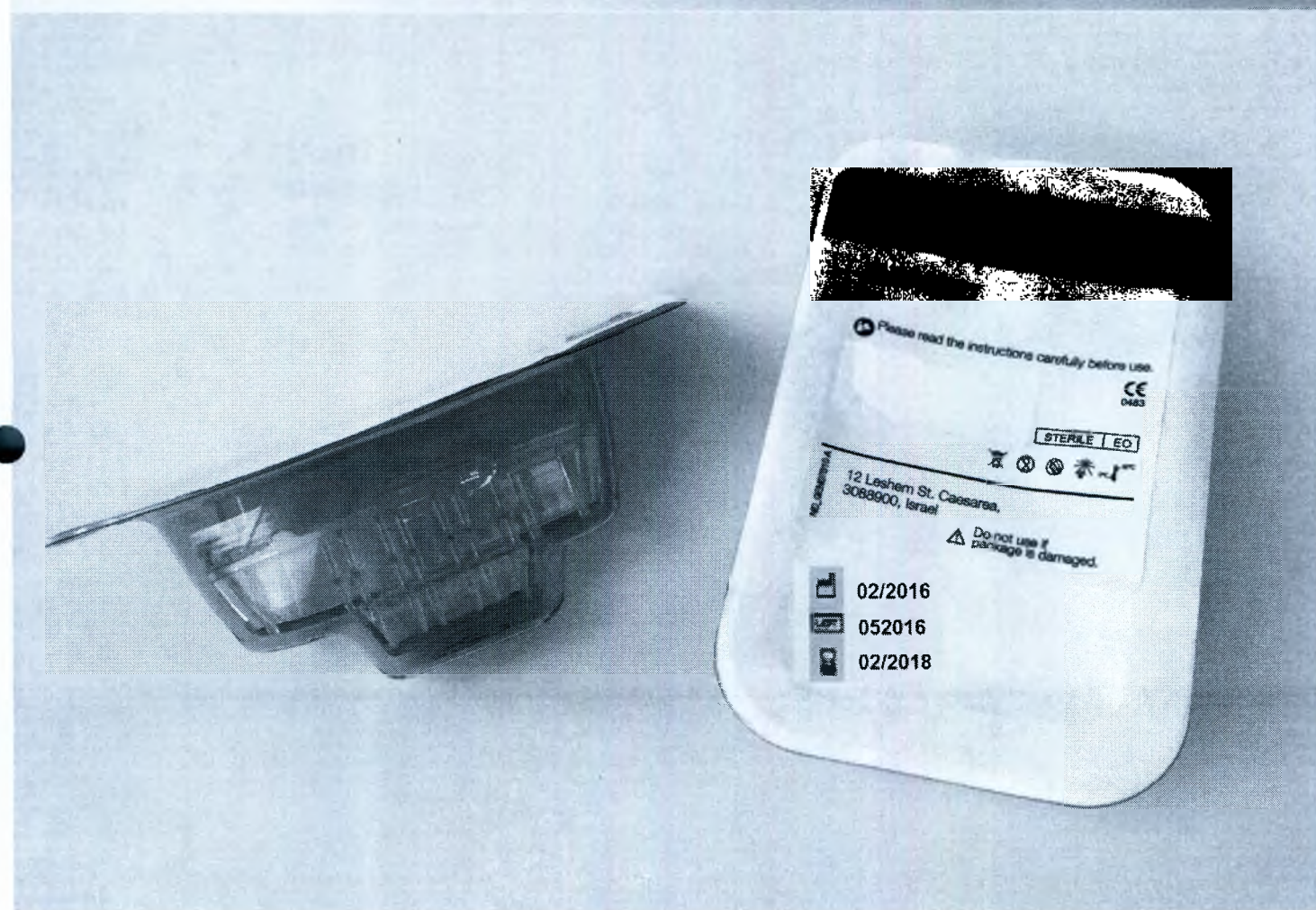1. Place the template on the intended treatment area.
2. Use a light, oil-free makeup base and a small sponge to spread the makeup through the template holes, creating a uniform matrix grid on the skin.
3. Position one corner of the lower edge of the FSR tip on one of the outside dots of the grid, so the surface of the tip covers the first area to be treated. Deliver pulse.
4. Lift the FSR Handpiece and position the same tip corner on the next grid dot so that the lower edge of the tip now covers the next area. Deliver pulse.
5. Repeat Step 4 for full and uniform coverage of treatment area.

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