3773

RuñA

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APPROVED BY

Technical Director &
Management Representative

Carlos Janer

27/01/ 2020

**The device is a radiofrequency hyperthermic
device, in options:**

Radiofrequency hyperthermic device

ELITE NS

Radiofrequency hyperthermic device

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2020

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Content

05	INTRODUCTION
05	- The radio-frequency hyperthermic device in versions
06	- Contact Information
06	- Appointment
07	- Operating principle
07	- Resistive mode (RES)
07	- Capacitive mode (CAP)
08	- Continuous modulation
09	KEY WARNINGS
11	INDICATIONS
11	CONTRAINDICATIONS
11	SIDE EFFECTS
12	PRODUCT DESCRIPTION
12	- Appearance of the product
26	- Product design
31	- Appointment of product elements
32	- Product controls
33	PRODUCT OPERATION
33	- Initial installation of the device
34	- Turning on
34	- shutdown
35	INTERFACE NAVIGATION
36	- Time Based Procedures
37	- Capacitive Time Based Procedures (CAP)
38	- Procedures based on time counting in resistive mode (RES)
39	- Procedures based on the amount of energy transferred
40	- Capacitive-based procedure (CAP) based on the amount of energy transferred
41	- Resistive mode (RES) based on the amount of energy transferred
42	- Use of procedures
42	- Saved procedures
43	—Creation of procedures
47	—Manage management
48	—Changing folders and programs
48	—Removing folders and programs
48	—Moving programs
49	—Copying programs
49	—Loading stored procedures

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50	INDIBA PREPARED PROCEDURES
51	SETTINGS MENU
52	MONITORING PROGRAM
52	HEATING A PLANE RETURN ELECTRODE
53	PRODUCT USE INFORMATION
54	ICONS
55	- Remote Control
56	PRODUCT ELEMENTS
56	- Types of electrodes
56	- Return electrodes
57	- Handle for CAP electrodes
57	- Handle for RES electrodes
59	- Manual
61	- Placing the active electrode and return electrode on the patient's body
62	- Guide to the use of active electrodes
62	- Capacitive electrodes (CAP)
63	- Resistive electrodes (RES)
65	- Cable Connection Guide
66	INDIBA ELECTRIC CONDUCTING CREAMS
67	PRODUCT CLEANING AND DISINFECTION
67	- Cleaning the device
67	- Cleaning and disinfection of active electrodes, return electrodes and electrode handles
67	- When carrying out gynecological procedures
68	REPAIR AND MAINTENANCE
69	MARKING AND PACKAGING
69	- marking
79	- Packing
86	TECHNICAL DETAILS
86	- Product classification
86	- Compliance with the standards of the Russian Federation
87	- Transportation, storage, operation
88	- Weight and overall characteristics
90	- Specifications
93	- Disposal
94	GUARANTEES
95	- Serial number
96	TROUBLESHOOTING AND TROUBLESHOOTING
98	- Error messages
99	GLOSSARY OF TECHNICAL TERMS
	EMC Information
	Appendix 1 Monitoring Program
	Annex 2 Protocols of Procedures

Introduction

HYPERTHERMAL RADIO FREQUENCY APPARATUS IN PERFORMANCE OPTIONS:

The device is radio frequency hyperthermic

In versions:

1. The device radio-frequency hyperthermic ELITE NS, as a part:

- The main unit of the device ELITE NS – 1 unit;
- Column – 1 unit;
- Shielded mains cable – 1 unit;
- Connection cable – 1 unit;
- Key – 6 unit;
- RES curved handle with cable (if required) – 1 unit;
- CAP curved handle with cable (if required) – 1 unit;
- CAP straight handle with cable (if required) – 1 unit;
- Remote control unit – (if required) – 1 unit;
- Cylindrical neutral electrode (if required) – 1 unit;
- Cable for cylindrical neutral electrode (if required) – 1 unit;
- Return plate (if required) – 1 unit;
- RES electrode 35 mm (if required) – 1 unit;
- RES electrode 50 mm (if required) – 1 unit;
- RES electrode 65 mm (if required) – 1 unit;
- RES electrode 90 mm (if required) – 1 unit;
- CAP electrode 30 mm (if required) – 1 unit;
- CAP electrode 40 mm (if required) – 1 unit;
- CAP electrode 55 mm (if required) – 1 unit;
- CAP electrode 65 mm (if required) – 1 unit;
- CAP electrode 80 mm (if required) – 1 unit;
- Manual – 1 unit;
- IDC NS letter download MS – 1 unit;

2. Welcome Pack, in composition:

- Proionic Anticellulite Cream – 1 unit;
- Proionic Face Cream – 1 unit;
- Proionic Care Cream – 2 units;
- Massager Electrode 65mm – 2 unit;
- USB – 1 unit;
- Bag for electrodes – 1 unit;

3. PACK GINECO-AESTHETICS, in composition:

- RES Longitudinal electrode (if required) – 1 unit;
- CAP Longitudinal electrode (if required) – 1 unit;
- RES electrode 35 mm (if required) – 1 unit;
- CAP electrode 30 mm (if required) – 1 unit;
- CAP electrode 40 mm (if required) – 1 unit;

Apparatus radio-frequency hyperthermic RGN, consisting of:

1. Radiofrequency hyperthermia device RGN, in composition:

- The main unit of the device RGN – 1 unit;
- Shielded mains cable – 1 unit;
- RES curved handle with cable (if required) – 1 unit;
- CAP curved handle with cable (if required) – 1 unit;
- Remote control unit (if required) – 1 unit;
- Return plate (if required) – 1 unit;
- RES electrode 35 mm (if required) – 1 unit;
- RES electrode 50 mm (if required) – 1 unit;
- CAP electrode 30 mm (if required) – 1 unit;
- CAP electrode 40 mm (if required) – 1 unit;
- CAP electrode 55 mm (if required) – 1 unit;
- Manual – 1 unit;
- IDC NS letter download MS – 1 unit;

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3. PACK GINECO-AESTHETICS, in composition:

- RES Longitudinal electrode (if required) – 1 unit;
- CAP Longitudinal electrode (if required) – 1 unit;
- RES electrode 35 mm (if required) – 1 unit;
- CAP electrode 30 mm (if required) – 1 unit;
- CAP electrode 40 mm (if required) – 1 unit;

For the sake of brevity, the radiofrequency device hyperthermic ELITE NS may be referred to as ELITE NS, and the radiofrequency device hyperthermic RGN - RGN.

CONTACT INFORMATION

Manufacturer: Indiba, S.A. («Indiba SA », Spain)
 Manufacturer Address: C/Moianés, 13, Polígono Industrial Can Casablanças, 08192 Sant Quirze del Vallés (Barcelona), Spain
 Phone: +34 93 265 55 22/+34 93 721 69 20;
 Fax: +34 93 576 30 40/+34 93 721 47 74
 Place of production: Indiba, S.A. ("Indiba SA"), Spain

Information about the authorized representative of the manufacturer in the Russian Federation:

BW Company Limited Liability Company (BW Company LLC)

Phone: +7 (499) 281-67-68

Email: info@beawire.com

Address: 129085, Moscow, Prospect Mira, 101, bldg. 1, pom. 17

PURPOSE

The product is intended for heating subcutaneous tissues due to radio frequency electromagnetic energy in order to reduce the appearance of cellulite, sagging skin, reduce wrinkles, and it is also used in postoperative tissue repair. The product can also be used to relieve pain, muscle spasms, including the treatment of chronic pelvic pain. The product consists of monopolar applicators connected to an energy generator. The product is not intended for the treatment of malignant tumors.

All-Russian product classifier by type of economic activity (OK 034-2014):

26.60.13.130

The class, depending on the degree of potential risk of medical use in accordance with the nomenclature classifier of medical devices, approved by order of the Ministry of Health of the Russian Federation dated June 06, 2012 No. 4n:

2a

Product classification in accordance with the nomenclature classifier of medical devices, approved by order of the Ministry of Health of the Russian Federation of June 06, 2012 No. 4n:

237180

Type of contact with the body: Radio-frequency hyperthermic apparatus - short-term contact (less than 24 hours) with intact skin; Longitudinal electrode RES, Longitudinal CAP electrode - short-term contact with intact mucous membranes.

Method of application: in accordance with the instruction manual.

Application conditions: medical institutions, used by qualified and specially trained specialists for adult patients.

OPERATING PRINCIPLE

In the line of devices manufactured by Indiba, S.A. ("Indiba S.A."), Spain (hereinafter referred to as INDIBA, Indiba) uses technology developed by INDIBA to non-invasively stimulate the healing and recovery mechanisms. The technology uses current circulation of a certain frequency (~ 448 kHz), also known as radio frequency current.

The device is a current generator with two operating modes: capacitive (CAP) and resistive (RES).

Current is supplied to the tissue using active electrodes (all electrodes except return electrodes) and a return electrode. To understand the principle of operation, it is necessary to imagine that the current generated in the product enters the active electrode, flows through the human tissues, reaches the return electrode, and then returns to the apparatus, closing the circuit.

The surface area of the active electrode is smaller than that of the return electrode, so the current only affects a certain surface area of the body. This means that the current density in the tissue under the active electrode is greater than the density in the tissue in contact with the return electrode.

In the operating modes described below (resistive and capacitive) at both electrodes (active and return), conductive creams made by INDIBA should always be used to ensure the necessary electrical contact between the electrodes and the tissue.

RESISTIVE MODE (RES)

In resistive mode, an active stainless steel electrode is used. This electrode is a good electrical conductor and acts directly on the treated tissue area (an electrically conductive cream is necessarily applied to the skin).

The current density in the tissue changes proportionally when flowing from the active electrode to the return electrode and reaches a maximum value near the active electrode. As a result of the current, the internal temperature of the tissue rises. The temperature reaches its maximum value near the active electrode, gradually decreasing as it moves away from it.

CAPACITY MODE (CAP)

In capacitive mode, a metal active electrode is used, coated with an insulating material, which acts directly on the treated tissue area (an electrically conductive cream is necessarily applied to the skin).

The electrode has a metal shaft "A" and a metal disk "B" coated with insulating material "C". The current leaving the apparatus passes through the electrode from the shaft "A" to the disk "B" and then is transmitted to the human tissue through the material "C".

As in resistive mode, current flows from the electrode through the tissue to the return electrode ("D"). In this case, the current increases the temperature of the tissue through which it passes. The increase in temperature and the depth of exposure depend on the voltage supplied by the device (depending on the power level selected by the user), the size of the electrode, the characteristics of the treated tissue area and the duration of the procedure.

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In capacitive mode, due to the insulating coating, heat is also generated in the lower part of the electrode, which makes it the preferred mode for processing subcutaneous areas.

In order for the current to pass through the insulation coating, the apparatus must generate a high output voltage.

Despite the fact that both methods increase the internal temperature of the tissue, higher temperatures and processing depths are achieved in the resistive mode.

CONTINUOUS MODULATION

The technology of frequency modulation was first used in surgical electroscalpels - it made it possible to control tissue coagulation and local bleeding.

Modulation is an additional functional mode of the device, which allows you to change the frequency amplitude, which has predefined characteristics. In this mode, the device automatically limits the output power value.

Modulation increases the speed and efficiency of tissue processing, allowing you to reduce the required power and improving the patient's subjective perception during treatment.

This mode can be turned on / off depending on the indications for treatment, the origin and duration of the problem, the location of the affected area, its degree and severity, the treatment regimen and, ultimately, in accordance with the therapeutic objectives.

Key warnings

IMPORTANT: Read the instruction manual carefully before using the machine.

Failure to follow the instructions described in the instruction manual may pose a risk to patients.

The device must be used by professional users who have completed a training course conducted by INDIBA.

WARNINGS:

- To avoid electric shock, the device should only be connected to a grounded electrical circuit.
- Make sure that the required voltage specifications are observed when connecting the device. The range of permissible supply voltages is indicated on the marking on the rear panel of the device.
- For the control panel, use only batteries with the characteristics specified in the instruction manual. Dispose of used batteries in accordance with current regulations.
- Leave a gap of at least 10 cm between the back of the unit and the wall to allow air to circulate through the fan. Make sure that the back of the device is not blocked - make sure that towels or clothes do not block the ventilation holes.
- Use only conductive creams made by INDIBA. The use of other electrically conductive creams may result in poor results and void the warranty on the device. Unapproved electrically conductive creams can cause burns and reduce the life of the electrodes.
- Use only cables supplied by the manufacturer. The use of other cables (for example, longer cables) can change the output values of the device or its resistance to other types of radiation.
- The device must be connected to a power source with an appropriate voltage. The range of permissible supply voltages is indicated on the marking on the rear panel of the device.
- To avoid electric shock, this unit should only be connected to a grounded power source.
- Never alter the design of the unit.
- Do not open or disassemble the unit. In the event of a malfunction: this unit does not contain parts that can be serviced by the user. Repairs to the device should only be carried out by qualified specialists certified by INDIBA.
- In the event of a malfunction or leakage of liquid from the liquid crystal display of the unit, avoid contact with the liquid.
- Do not use the product near short-wave therapy devices, as their radiation may affect the product.
- This unit is not intended for use in flammable environments.
- This device should not be used on patients at the same time as any other product.
- Do not use the device on patients with electronic implants (such as a pacemaker) without a physical examination.
- There should be no rings, earrings, piercings, or other metal objects on the area of the patient to be treated. Avoid contact with conductive parts that may create unwanted paths for radio frequency current. When using the machine, the operator must remove the rings from all fingers.
- The operator must not touch the electrodes with his hands. Always hold the electrode by the electrode handle.
- If the unit is operating, press the pause button before inserting, removing or replacing the electrode.
- Never use capacitive electrodes with signs of wear, such as small holes, cracks, or chips.
- Never use the unit if the return electrode has not been properly secured. Failure to do so may result in burns to the patient.
- Do not use electrodes if there are signs of damage or wear on the cables of the electrode handle.
- Stop the procedure immediately if the patient feels that the electrodes are too hot.
- Before starting the procedure, clean the treated area.
- Remove any products or substances that may affect the procedure from the skin.

- Do not apply excessive force to the shafts of the electrodes when using them; improper handling can lead to their displacement or cracking on the surface of the electrode.
- CAP electrodes are coated with polyamide plastic. Do not expose them to temperatures above 90°C.
- Avoid impact of electrodes, do not scratch them with sharp or pointed objects. If the coating is worn or breaks off the metal part, do not use the CAP electrode - replace it.
- Hypotensive patients undergoing treatment for the cervical spine or any other part of the spine or joints of the limbs may experience a decrease in blood pressure. In this case, treatment should be stopped until the blood pressure level returns to normal.
- Before performing facial treatments in resistive mode, the false jaw should be removed.
- Direct exposure to the testes should only be performed at subthermal temperature.
- Patients taking oral anticoagulants should consult a doctor before starting treatment with this unit.
- The device is intended only for the treatment of adult patients.
- In elderly patients, skin lesions are more likely to occur, therefore treatment should be carried out with extreme care to avoid injuries or skin burns.
- The device should be used with caution in patients with known allergic reactions to nickel and chromium. Stop treatment if an allergic reaction is observed in areas in contact with the return flat electrode, return cylindrical electrode or resistive electrode, or electrically conductive creams.
- The use of resistive treatment regimens in areas with internal infection or encapsulation (abscesses) is not recommended. There is a risk of infection. In this case, you must consult your doctor before starting treatment with this unit.
- When performing gynecological procedures, a specialist is advised to use nitrile or latex disposable gloves (not included in the delivery).

Indications

Skin Tightening:

- Reduction of wrinkles.
- Strengthening the neck and décolleté.
- Reduction of double chin.

Cellulite reduction

- Gynoid lipodystrophy

Body remodeling

- Non-invasive figure correction
- Removal or reduction of excess skin and fat
- Body contouring

Before and after surgery

- Preoperative treatment
- Tissue healing
- Wound healing
- Getting rid of scars
- Burn treatment

Vaginal rejuvenation and related clinical problems

- Correction of vulvovaginal weakness
- Treatment of stress urinary incontinence
- Genitourinary syndrome in menopause
- Sexual dysfunction.

The product can also be used to relieve pain, muscle spasms, including the treatment of chronic pelvic pain. Any treatment, whether monotherapy or a combination of different types of treatment, should be carried out only as prescribed by the doctor.

Contraindications

- The presence of pacemakers or other electronic implants.
- Pregnancy.
- Skin lesions (open wounds or recent burns).
- Thrombophlebitis.
- Lack of sensitivity (congenital insensitivity to pain, nerve damage, paraplegia or pharmacological treatment that reduces sensitivity to pain or heat)
- The product is not intended for the treatment of malignant tumors.

Side effects

After the procedure, patients usually feel warm in the treated area, which may cause some discomfort. In rare cases, the use of electrically conductive creams can cause slight skin irritation.

Product Description

APPEARANCE OF THE PRODUCT

Fig. 2 General view of the product in the embodiment. Radio-frequency hyperthermic device ELITE NS (the main unit is located on the trolley-rack, the drawers are extended, the handles are curved for RES and CAP electrodes and a return flat electrode)

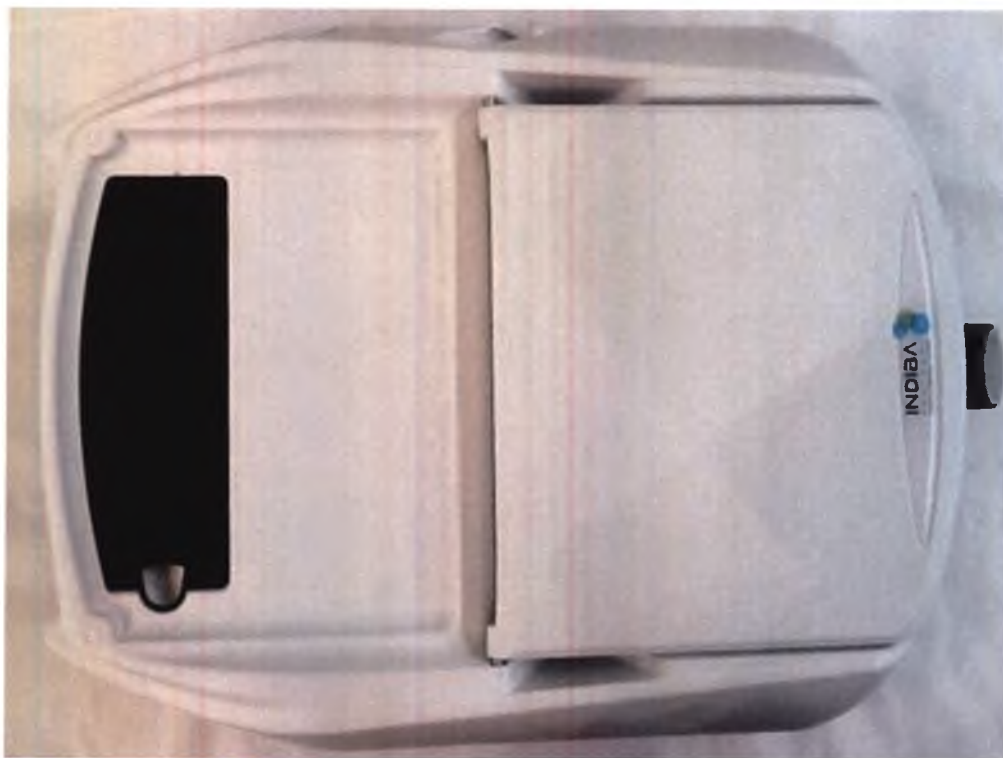
Fig. 3 Appearance of the main unit in the ELITE NS and RGN version (top view with the lid closed)

Fig. 4 Appearance of the main unit in the ELITE NS and RGN version (top view with open lid)

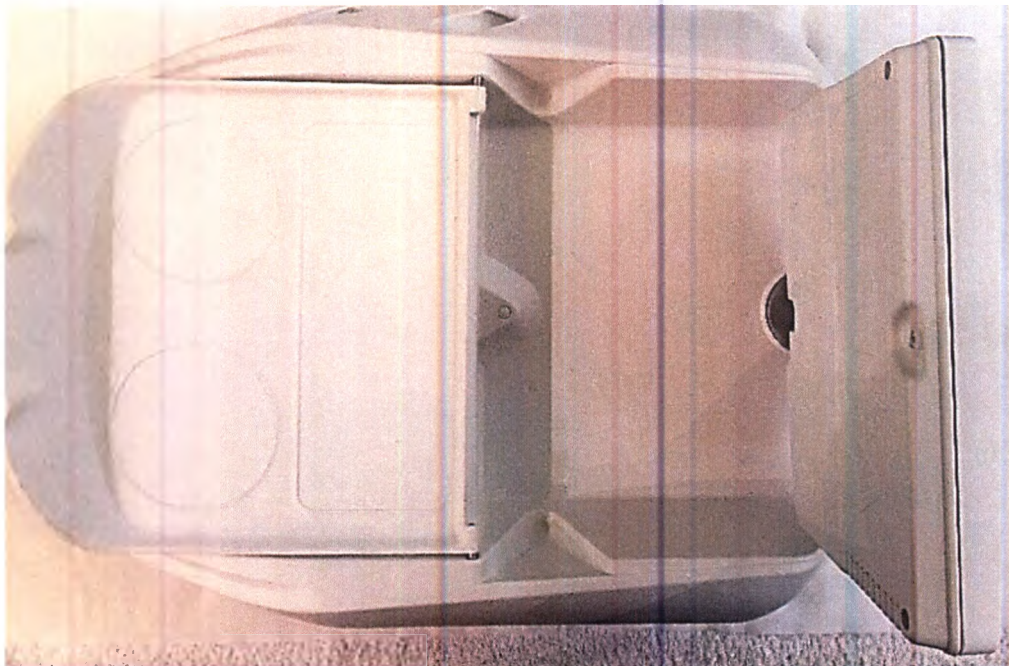
Fig. 5 Appearance of the main unit in the ELITE NS and RGN version (top view with the cover open and the monitor raised)

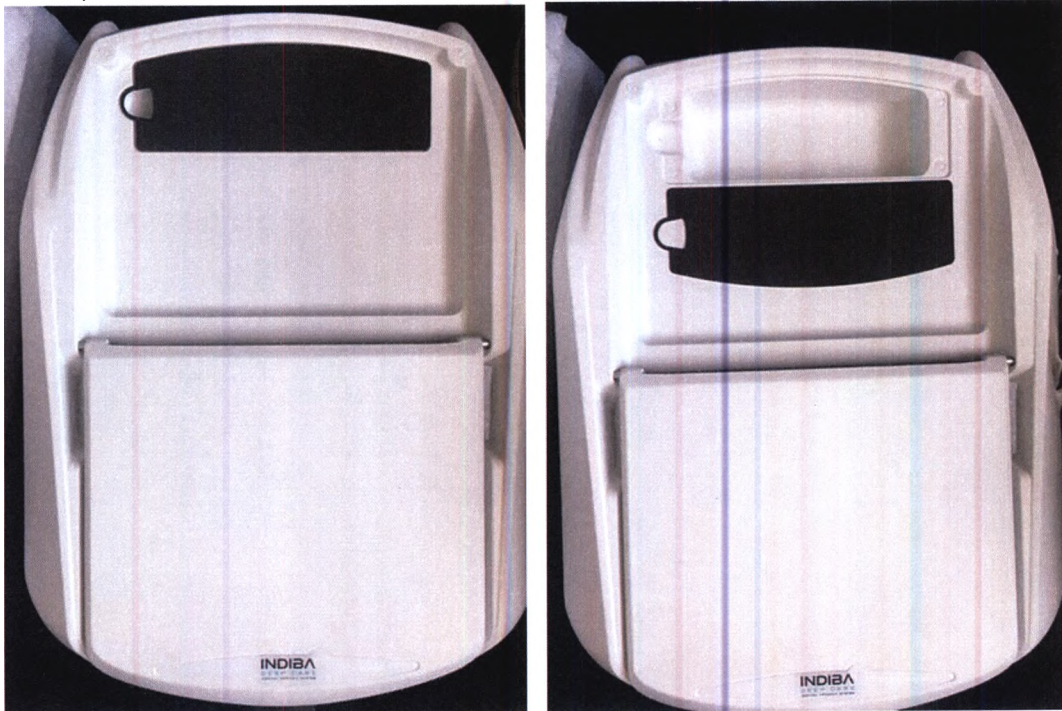
Fig. 6 Product appearance with a closed (left) and open (right) compartment for the control panel. The cover for the control panel is mounted on the housing using built-in magnets.

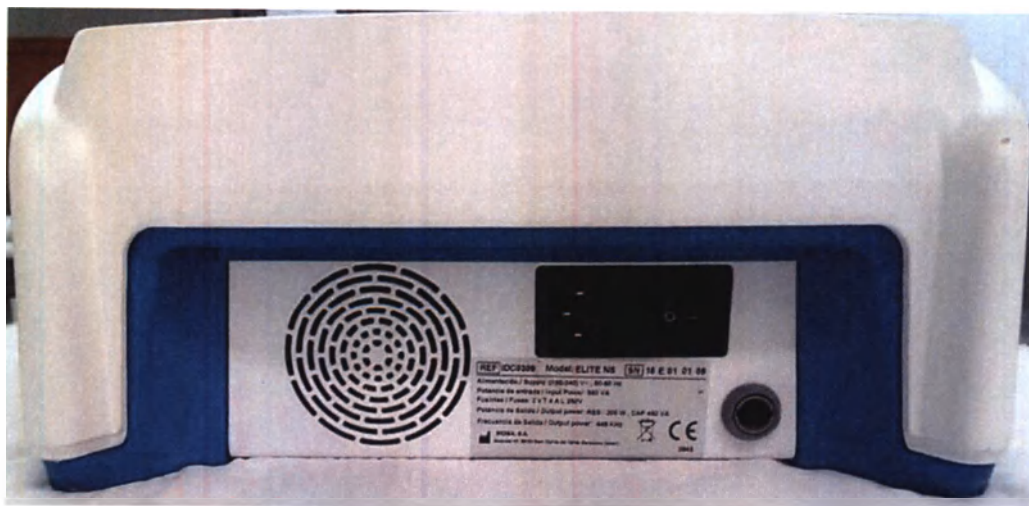
Fig. 7 Appearance of the main unit in the ELITE NS and RGN version (rear view, marking is presented for the ELITE NS version)

Fig. 8 Appearance of the main unit in the ELITE NS and RGN version (side view)

Fig. 9 Socket for connecting return flat and cylindrical electrodes. In order to remove the connector for connecting to the device from the device itself, you must press the top "PUSH" button

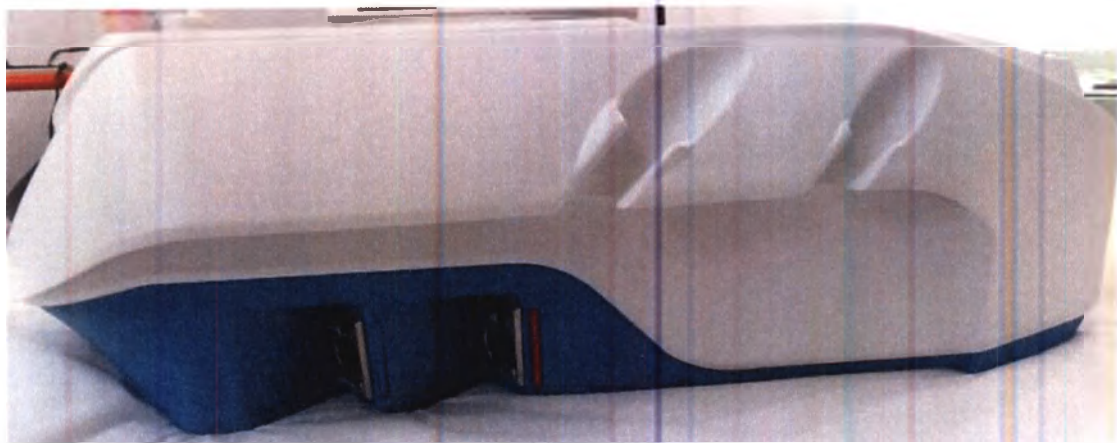
Fig. 10 Appearance of the main unit in the ELITE NS and RGN version (side view)

Fig. 11 Connector for connecting handles for electrodes. Next to each connector is a specific color: "blue" - this connector is for the handle curved for the RES electrode (there is a blue bezel on the handle as a hint to the operator), "red" - this connector is for the handle curved and straight for electrodes CAP (there is a red bezel on the handle as a hint to the operator). In order to remove the connector for connecting to the device from the device itself, you must press the top "PUSH" button.

Fig. 12 Appearance of the main unit in the ELITE NS and RGN version (front view)

Fig. 13 Appearance of the main unit in the ELITE NS and RGN version (bottom view)

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Fig. 15 Shielded power cable

Fig. 16 Bag for electrodes

Fig. 17 Open bag for electrodes. The bag for electrodes opens and closes with a sewn-in zipper. Foamed polyethylene is inserted into the bag to place the electrodes (in the center of each black circle from the pls there is a special hole for the electrode shaft).

Fig. 18 Connection cable

Fig. 19 Curved handle for CAP electrode

Fig. 20 Curved handle for RES electrode

Fig. 21 Straight handle for CAP electrode

Fig. 22 flat return electrode

Fig. 23 Cylindrical return electrode

Fig. 24 Cable return electrode cylindrical

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Sección 2ª

[Signature]

Fig. 25 RES electrode (left) and CAP (right) (also referred to as resistive and capacitive electrodes, respectively)

Fig. 26 Conductive Cream Care Cream (left), Anticellulite Cream (center), Face Cream (right) (in a plastic protective bag)

Fig. 27 USB stick

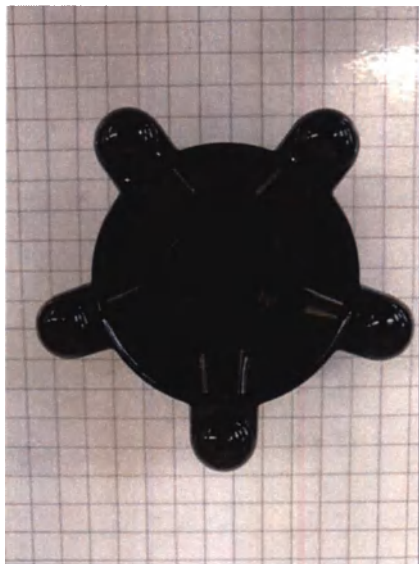
Fig. 28 Massage nozzle

Fig. 29 Trolley (front view) with keys inserted.

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Fig. 30 Rack-trolley (side view with connected cable connecting)

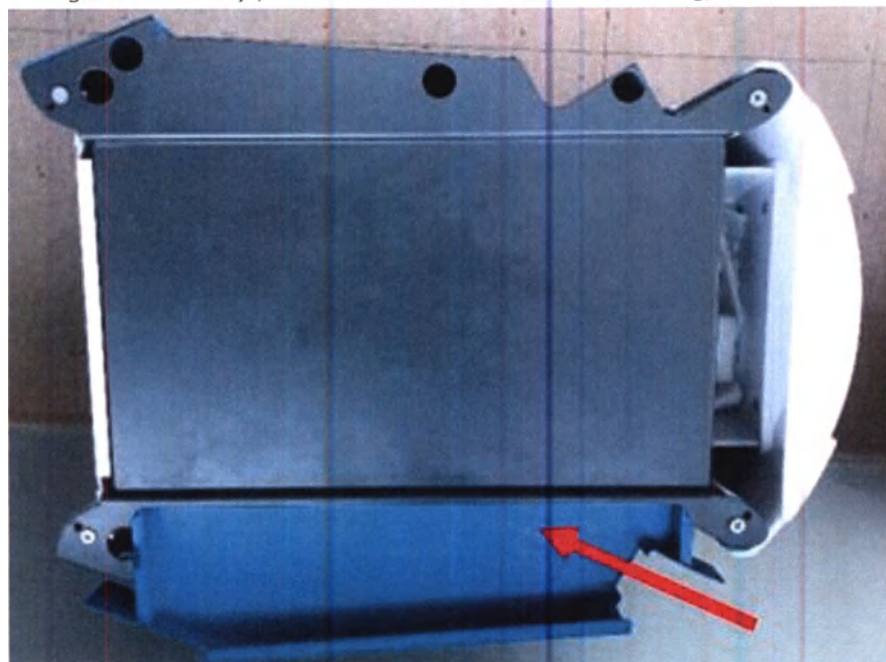
Fig. 31 Rack-cart (top view). The arrow indicates the retractable compartment to accommodate the return flat electrode

Fig. 32 Voltage switching at the rear of the trolley rack and connector for connecting the main unit with the trolley rack

Fig. 33 Longitudinal CAP electrode

Fig. 34 RES electrode longitudinal

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
PRODUCT DESIGN


Front View

1. Touch screen (display)
2. Handle connector for RES electrodes
3. Handle Connector for CAP Electrodes
4. Return electrode connector
5. Power and trigger settings knob / Stopping procedure (also called encoder)

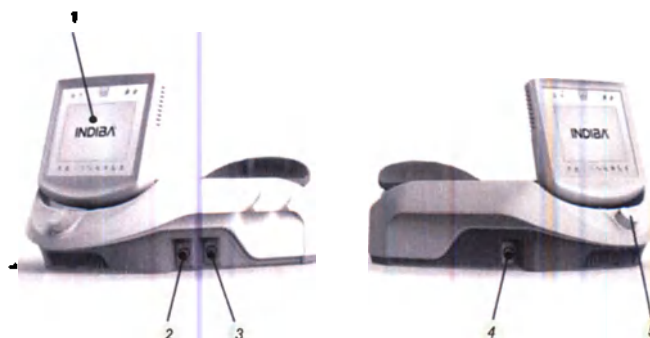Rear View


6. Power switch
7. Fuse compartment
8. Shielded power cable connector
9. Fan
10. Marking
11. Jack for connecting the cable of the connecting rack trolley

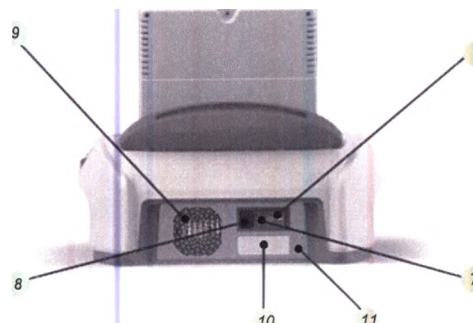Remote Control

12. Power down button
13. Power increase button
14. Button to start / stop the procedure

The status indicator is marked with an asterisk: on - the device has detected contact with the patient (the procedure is in progress), blinking - the device is waiting for contact with the patient, off - the button for starting the procedure has not been pressed.

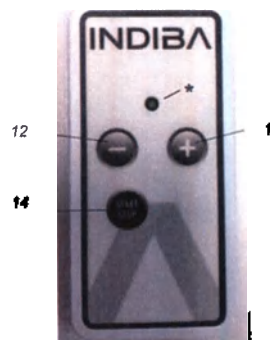* Comes with a battery for the remote control.

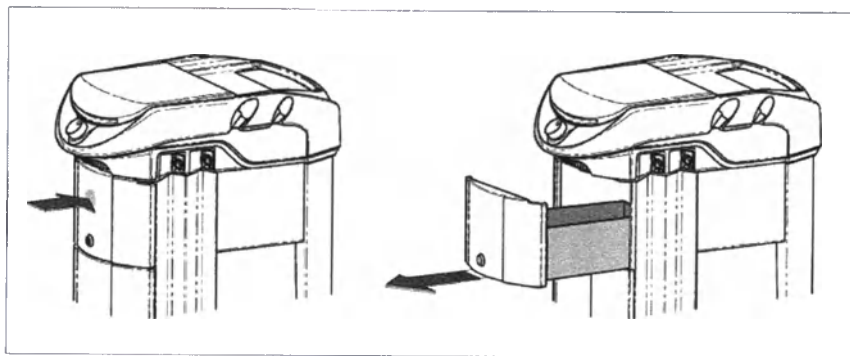
Fig. 35 Opening the drawer racks

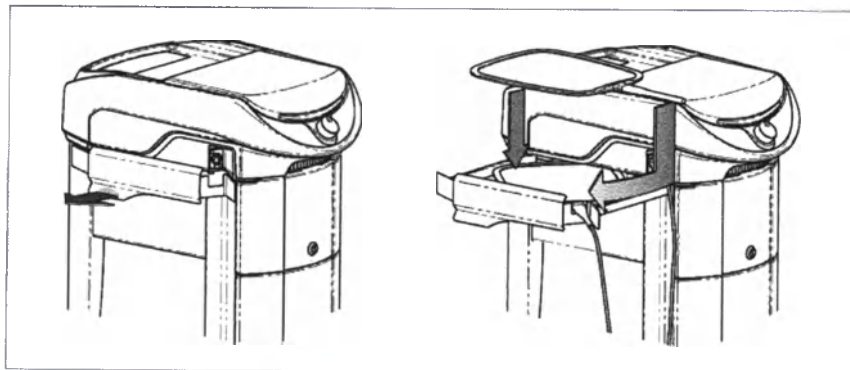
Fig. 36 Opening the drawer for the electrode return flat

Fig. 37 The locking mechanism of the wheels. It is necessary to press the foot on the elements marked by arrows to block the movement of the rack-cart

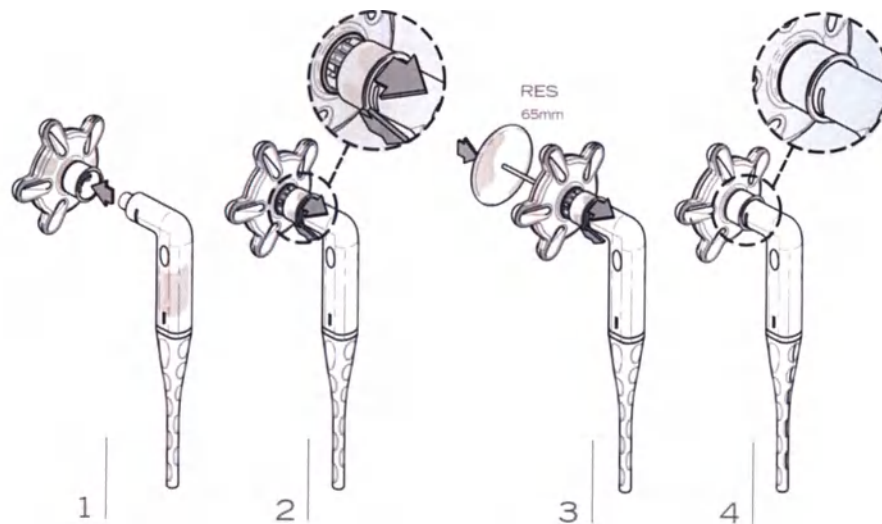
Fig. 38 Installation instructions for the massage nozzle. The massage nozzle can only be used with a 65 mm diameter RES electrode.

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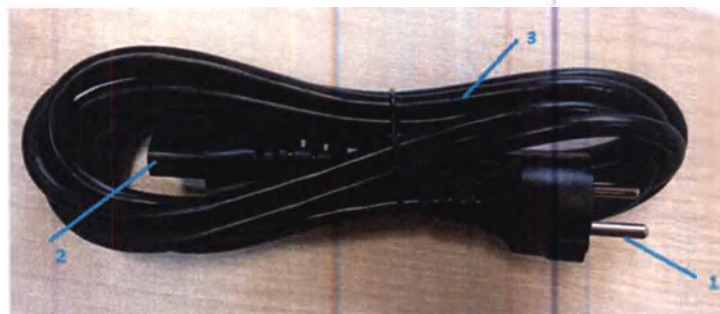
Fig. 39 Shielded Power Cable Design

- 1 - Power plug for connecting to the power supply
- 2 - Connector for connecting to the device
- 3 - Shielded power cable wire

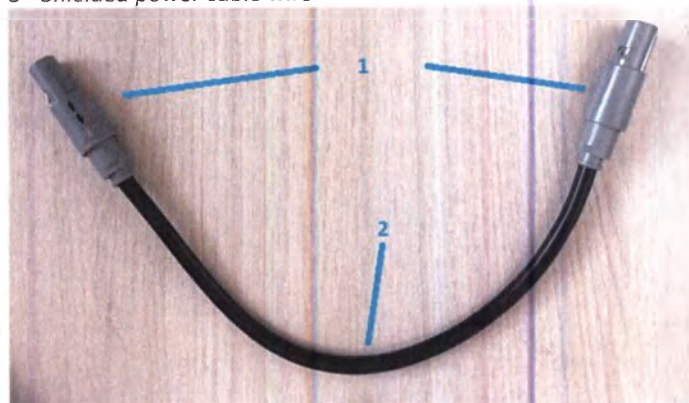
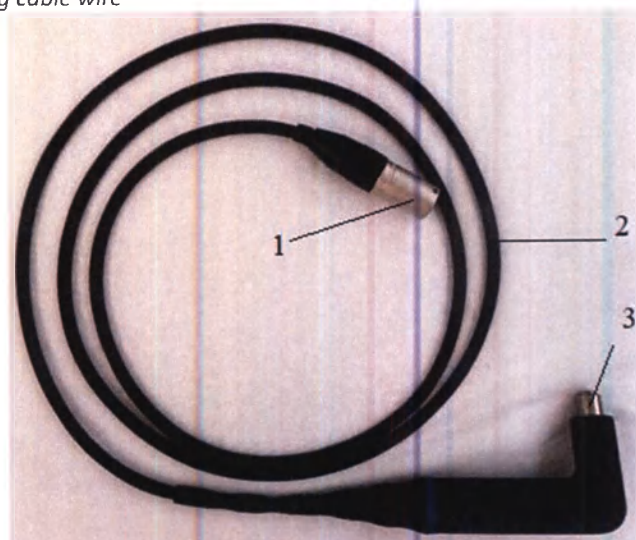
Fig. 40 Cable design

- 1- Connectors for connecting to the trolley-rack and the main unit (have the same interface)
- 2 - Connecting cable wire

*Fig. 41 Curved handle design for RES electrode
(other handles have a similar design)*

- 1 - Connector for connecting to the device
- 2 - The wire of the handle bent for the electrode RES
- 3 - Connector for connecting the RES electrode

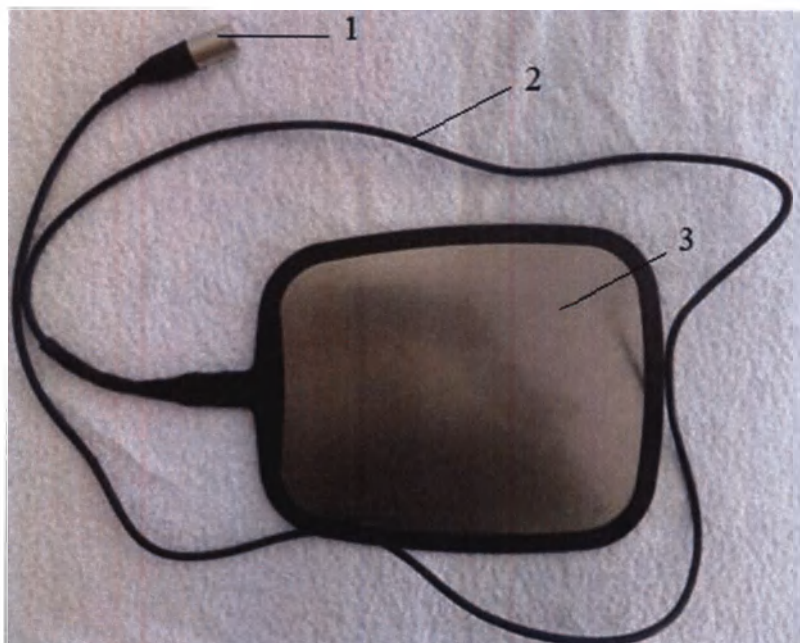
Fig. 42 Flat electrode return design
 1 - Connector for connecting to the device
 2 - wire return flat electrode
 3 - Return flat electrode

Fig. 43 The design of the cable return electrode cylindrical
 1 - Connector for connecting to the electrode return cylindrical
 2 - Connector for connecting to the device
 3 - Wire cable electrode return cylindrical

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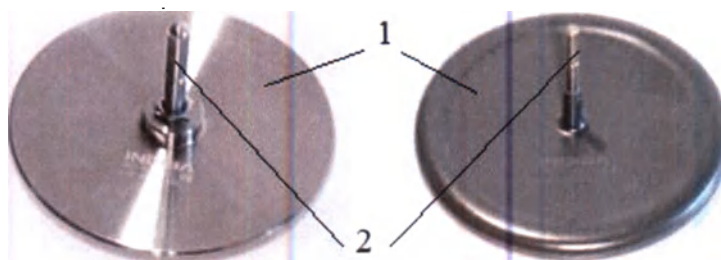
Fig. 44 Electrode design RES (left) and CAP (right)

1. The electrode
2. Electrode shaft

APPOINTMENT OF PRODUCT ELEMENTS


Block main ELITE NS, Main RGN Block	Generates a high-frequency electric current, controls the operation of the device.
Trolley	Serves as a stand for the main unit and the method of its movement using the wheels of the trolley rack. Its drawers can be a storage place for product items.
Shielded power cable	Used to connect the device to the power supply.
Connecting cable	Serves to transfer energy from the main unit to the trolley rack. Since the return flat electrode is in direct contact with the patient, it is important that the patient is comfortable touching him. Inside the trolley rack there is a heater, and the return flat electrode placed in a special compartment is heated to a comfortable temperature.
Cylindrical Return Electrode Cable	Using this cable, a return cylindrical electrode can be connected to the device.
Electrode bag	In it you can store and transfer electrodes.
Key	Trolley racks can be locked with keys.
Handle (all versions)	Used to connect electrodes
Electrode RES, CAP (all versions)	Serves for supplying radio frequency current to the patient's body.
Return electrode (all versions)	It serves to close the circuit between the patient and the main unit.
Пульт управления	Serves for remote control of the device
Conductive cream (all versions)	Serves to reduce the electrode-skin contact resistance during the procedure, also minimizes friction between the skin and the electrode.
Massage nozzle	Attaches to the electrode to create an additional massage effect during the procedure.
USB stick	The USB drive contains photos of the device, photos of the procedures, advertising photos and videos, as well as an instruction manual.

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Handwritten signature

PRODUCT CONTROLS

- 15. Return button to the main screen
- 16. Button to return to the previous screen
- 17. Access button in the settings menu
- 18. Button on / off monitoring program
- 19. On / off button for heating the return flat electron *
- 20. Product operation information display button

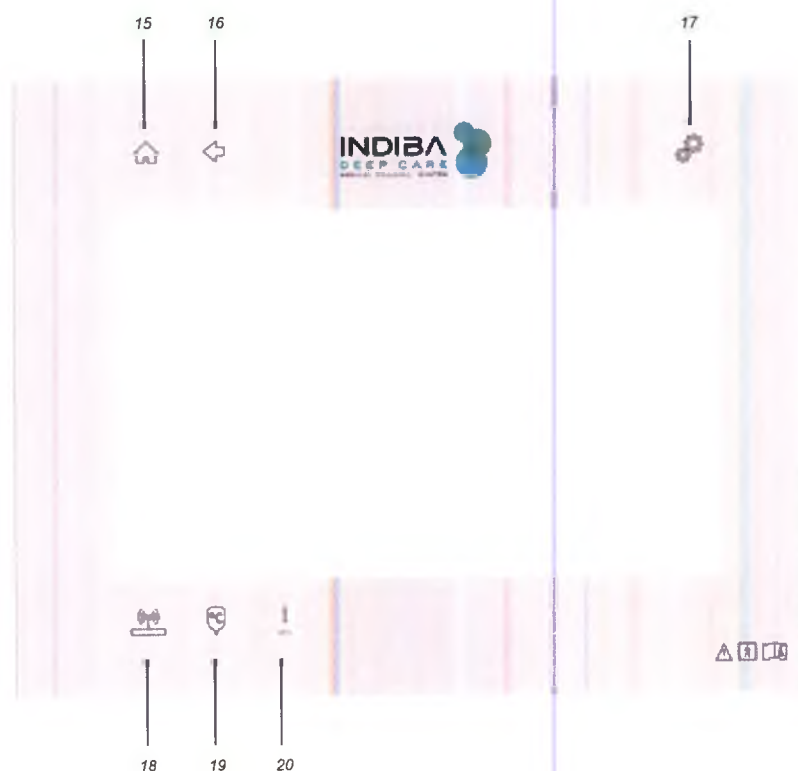The presence of these icons depends on the status of the display.


(*) This function is only available if the product is mounted on a trolley rack and connected to it using a connecting cable.

Product Operation


PRIMARY INSTALLATION

The initial installation of the product can be carried out by the end user.

In order to install the main ELITE NS unit on a trolley rack, it is necessary to adhere to the installation diagram shown below. On the block there are mainly special protruding elements, which, when installed correctly, should completely coincide with the holes on the top panel of the trolley-rack. A correctly installed main unit must not wobble on a cart stand.

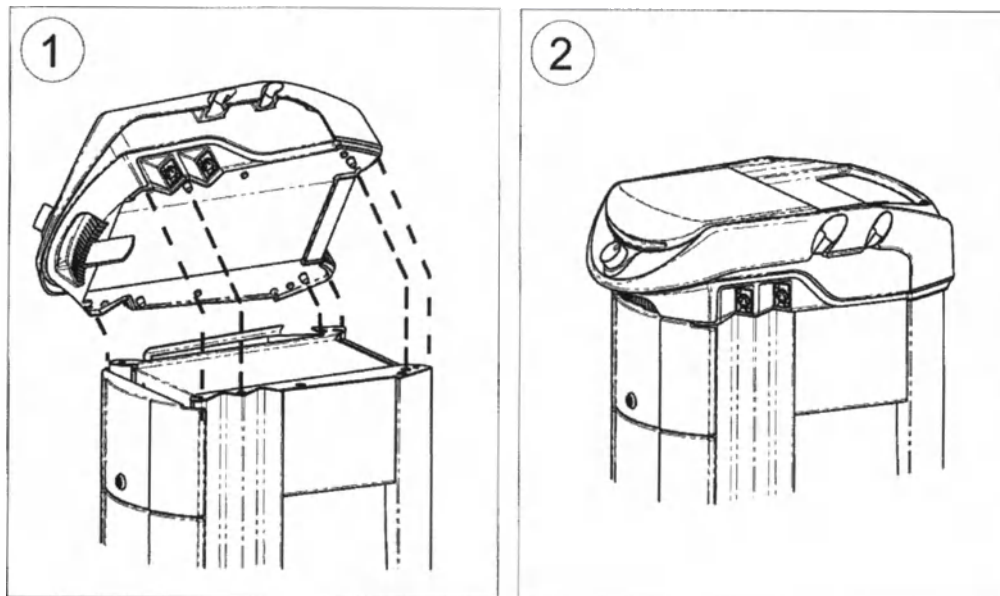
Before turning on the device, select the appropriate mains voltage.

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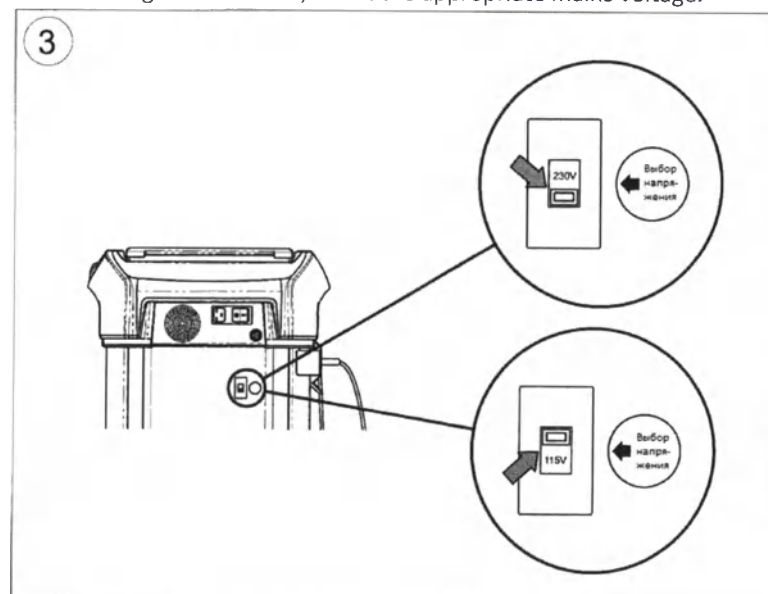Connect the connecting cable to the connector of the main unit to enable heating of the surface intended for placement of the return flat electrode.

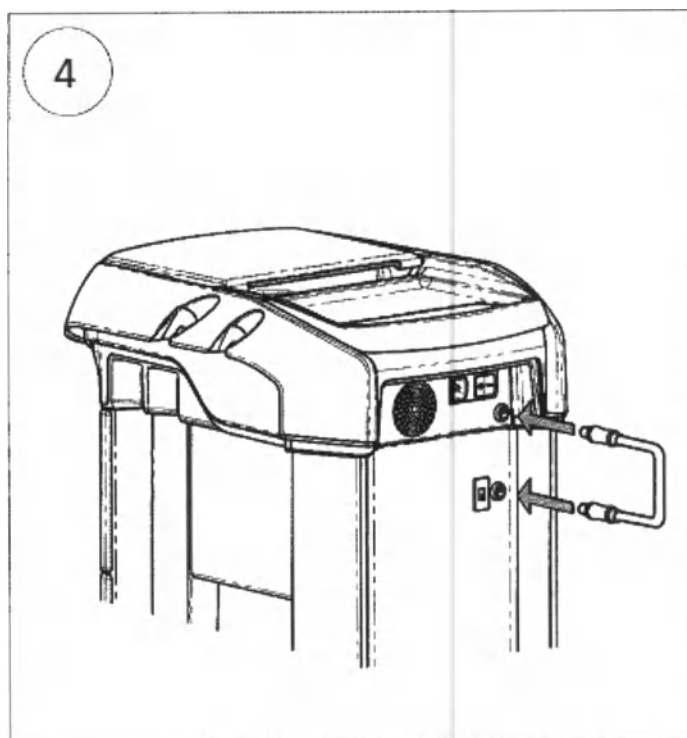TURNING ON

1. Connect the shielded power cable to the back of the unit (8).
2. Plug the shielded power cable into a power outlet.
3. Turn on the machine by pressing the power switch on the rear panel (6).

TURNING OFF

If the machine is not in use, press the power switch on the rear panel (6) and set it to the "off" position. It is recommended to turn off the product if it is not planned to use it for several hours.

Interface navigation

When turned on, the start screen should appear for a few seconds on the screen of the device (the device is loaded and checked at the same time).

At the end of the launch, the main screen should appear:

- 21. Button to start the procedure based on the countdown
- 22. Start button for a procedure based on the amount of energy transferred
- 23. Button to go to saved procedures

Buttons 21, 22 and 23 can be selected using the encoder (5).

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
[Handwritten signature]

TIME-BASED PROCEDURES

When you press the start button of the procedure based on the countdown (21), the device displays the following screen:

24. Button for increasing the duration of the procedure
25. Button to reduce the duration of the procedure
26. On / off button for monitoring the limit value of the output power
27. On / off button for modulating the output signal
28. Mode select button (CAPACITY (CAP) / RESISTIVE (RES))
29. Button to complete the procedure
30. "PAUSE" button, which pauses the procedure
31. The "START" button, which starts the procedure

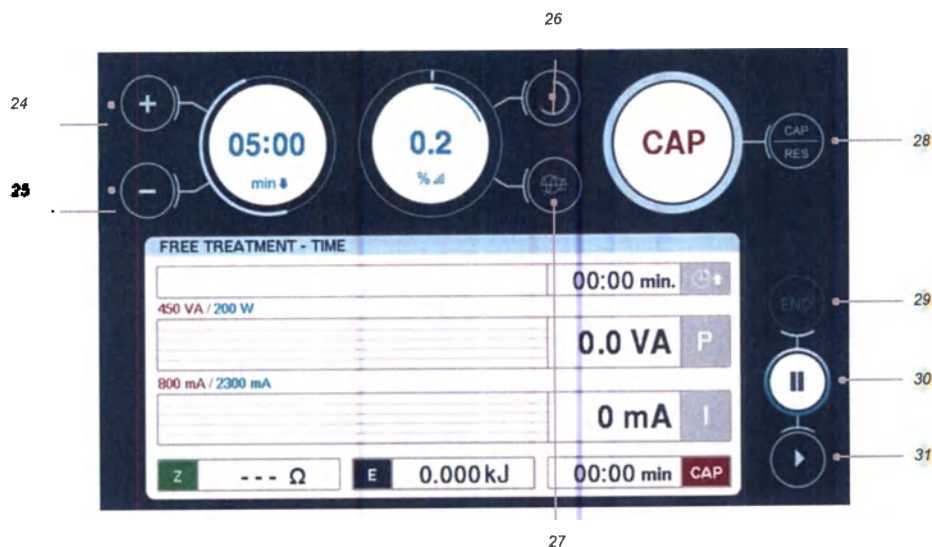Before starting the procedure, you can set the following values:

- Procedure time: to do this, press the increase button (24) and decrease the duration of the procedure (25). The minimum time is 1 minute and the maximum is 99 minutes. The duration of the procedure can be increased or decreased in increments of one minute.
- The output power can be changed using the encoder (5) or buttons (12) and (13) on the control panel. The output power can be changed from 0 to 100%, in each range of values its own increment step:
 - Output power up to 5% - step 0.2% for the ELITE NS model and 1% for the RGN model.
 - Output power from 5% to 30%, step 1%.
 - Output power from 30% to 100% - step 2%.
- Selecting an operating mode using the mode selection button (28).
- The control of the output power limit value (26) can be turned on / off in the "Settings" menu (17). This button (26) is used to enable or disable this option. You can change the limit value of the output power without going to the "Settings" menu by pressing and holding the button (26).
- Button (27) allows you to enable / disable the modulation of the output signal. When modulating, the device automatically limits the output power.

The following values can be set during the procedure.

- Output power - use the encoder (5) or buttons (12) and (13) on the control panel.
- Procedure time - you can increase (24) and decrease (25) the procedure time.

Capacitive Time Based Procedure (CAP)

The device displays the following operating parameters:

32. Time (in minutes) set by the user. After the procedure time has elapsed (the set time is displayed by the timer), the device emits a sound signal that marks the end of the procedure and stops supplying output power. A message appears on the display that the procedure has been completed.
33. The indicator of the output power level in%.
34. Operation mode indicator (CAP).
35. Time (in minutes). The cumulative time from the start of the procedure to the current moment.
36. Transmitted apparent power: current value of apparent power supplied in VA (volt-amperes). This value is the sum of resistive power (active or true) and reactive power (isolating). The main graph shows the change in this value over the past procedure time.
37. Output current (I): current value of the current flowing through the electrode and tissue, measured in milliamperes (mA). In the "Settings" menu (17) you can select the value (Ω or mA), the change of which you want to show on the graph.
38. Impedance (Z) measured in Ohms (Ω): current resistance value of body tissue and electrode. If three dashes ("---") appear, this value is too large and therefore goes beyond the scale, or the output power is too low to correctly perform this measurement. The product can be configured so that the second graph shows the change in this value during the procedure.
39. Energy (E), measured in kilojoules (kJ): the cumulative value of the power supplied during this step of the procedure.
40. Time (in minutes). The cumulative value of the duration of the procedure.


Time Based Resistive Procedure (RES)

The product displays the following operating parameters:

32. Time (in minutes) set by the user. After the procedure time has elapsed (the set time is displayed by the timer), the device emits a sound signal that marks the end of the procedure and stops supplying output power. A message appears on the display that the procedure has been completed.
33. The indicator of the output power level in%
34. Operating mode indicator (RES).
35. Time (in minutes). The cumulative time from the start of the procedure to the current moment.
36. Transmitted apparent power: current value of apparent power supplied in watts (W). The main graph shows the change in this value over the past procedure time.
37. Output current (I): current value of the current flowing through the electrode and tissue, measured in milliamperes (mA). In the "Settings" menu (17) you can select the value (Ω or mA), the change of which you want to show on the graph.


- 38 Electrical resistance (R), measured in Ohms (Ω): current value of tissue resistance. If three dashes ("---") appear, this value is too large and therefore goes beyond the scale, or the output power is too low to correctly perform this measurement. The product can be configured so that the second graph shows the change in this value during the procedure.
- 39 Energy (E), measured in kilojoules (kJ): the cumulative amount of energy supplied during this step of the procedure.
40. Time (in minutes). The cumulative value of the duration of the procedure.

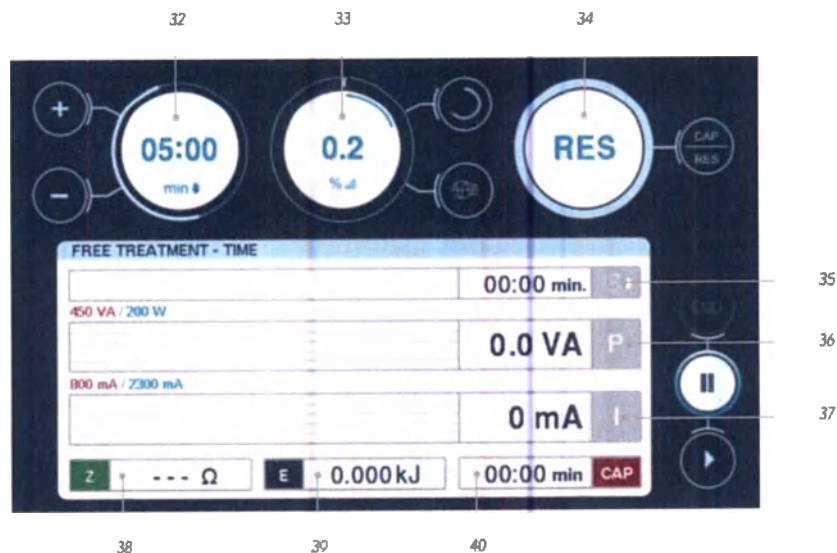PROCEDURES BASED ON THE AMOUNT OF TRANSFERRED ENERGY

When you press the start button of the procedure based on the amount of energy transferred (22), the device displays the following screen:

24. Button to increase the amount of energy transmitted.
25. Button to reduce the amount of energy transmitted.
- 26 On / off button for monitoring the limit value of the output power.
- 27 Button to enable / disable modulation of the output signal.
28. Mode select button (CAPACITY (CAP) / RESISTIVE (RES)).
29. Button to complete the procedure.
30. The "PAUSE" button, which pauses the procedure.
31. The "START" button, which starts the procedure.

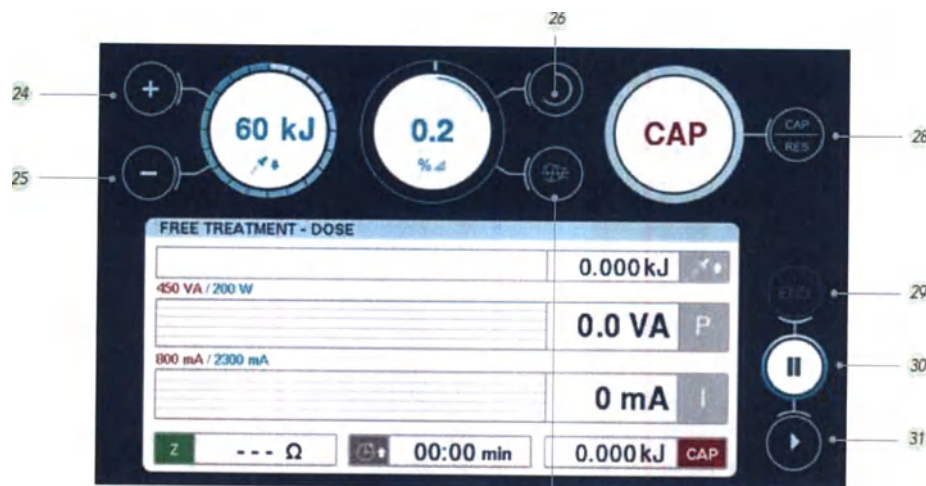Before starting the procedure, you can set the following values:

- The amount of energy transferred can be increased or decreased using the increase (24) and decrease (25) buttons. The minimum amount is 1 kJ, and the maximum is 999 kJ. The amount of energy transferred can be increased or decreased in increments of 1 kJ. You can also increase or decrease the value by 10 kJ by pressing and holding the button for the amount of energy transferred.
- The output power can be changed using the encoder (5) or buttons (12) and (13) on the control panel. The output power can be changed from 0 to 100%, in each range of values its own increment step:
 - Output power up to 5% - step 0.2% for the ELITE NS model and 1% for the RGN model.
 - Output power from 5% to 30%, step 1%.
 - Output power from 30% to 100% - step 2%.
- Selecting an operating mode using the CAP / RES button (28).
- The control of the output power limit value (26) can be turned on / off in the "Settings" menu (17). This button (26) is used to enable or disable this option. You can change the limit value of the output power without going to the "Settings" menu by pressing and holding the button (26).
- Button (27) allows you to enable / disable the modulation of the output signal. When modulating, the device automatically limits the output power.

During the procedure, the following values can be set:


- Output power - use the encoder (5) or buttons (12) and (13) on the control panel.
- Procedure time - you can increase (24) and decrease (25) the procedure time.

Capacitive Capacitive Transfer Procedure (CAP)

The product displays the following operating parameters:

32. The amount of energy transmitted (kJ) as set by the user. After reaching a predetermined amount of transmitted energy (the set value displays the counter of transmitted energy), the device emits a sound signal that marks the end of the procedure and stops supplying output power. A message appears on the display that the procedure has been completed.
33. The indicator of the output power level in%.
34. Operation mode indicator (CAP).
35. Energy (E), measured in kilojoules (kJ): the cumulative value of the energy transferred from the beginning to the current moment.
36. Total power: current value of the total power supplied in VA (volt-amperes). This value is the sum of resistive power (active or true) and reactive power (isolating). The main graph shows the change in this value over the past procedure time.
37. Output current (I): current value of the current flowing through the electrode and tissue, measured in amperes (A). In the "Settings" menu (17) you can select the value (Ω or mA), the change of which you want to show on the graph.


38. Impedance (Z) measured in ohms (Ω): current resistance value of body tissue and electrode. If three dashes ("---") appear, this value is too large and therefore goes beyond the scale, or the output power is too low to correctly perform this measurement. The product can be configured so that the second graph shows the change in this value during the procedure.
39. Time (in minutes). The cumulative time from the start of the procedure to the current moment.
40. Energy (E), measured in kilojoules (kJ): the cumulative amount of energy supplied during this step of the procedure

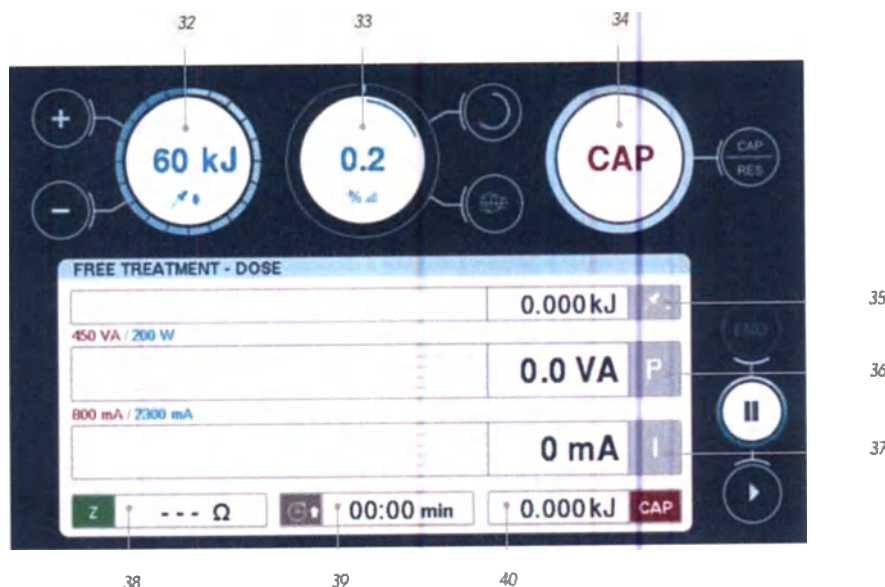Energy Resistant Transmission Procedure (RES)

The device displays the following operating parameters:

32. The amount of energy transmitted (kJ) as set by the user. After reaching a predetermined amount of transmitted energy (the set value displays the counter of transmitted energy), the device emits a sound signal that marks the end of the procedure and stops supplying output power. A message appears on the display that the procedure has been completed.
33. The indicator of the output power level in%.
34. Operating mode indicator (RES).
35. Energy (E), measured in kilojoules (kJ): the cumulative amount of energy supplied from the beginning of the procedure to the current moment.
36. Gross Power: The current value of gross power supplied in watts (W). The main graph shows the change in this value over the past procedure time.
37. Output current (I): current value of the current flowing through the electrode and tissue, measured in amperes (A). In the "Settings" menu (17) you can select the value (Ω or mA), the change of which you want to show on the graph.
38. Electrical resistance (R), measured in Ohms (Ω): current value of tissue tissue resistance. If three dashes ("---") appear, this value is too large and therefore goes beyond the scale, or the output power is too low to correctly perform this measurement. The product can be configured so that the second graph shows the change in this value during the procedure.
39. Time (in minutes). The cumulative time from the start of the procedure to the current moment.
40. Energy (E) measured in kilojoules (kJ): the cumulative value of the power supplied during this step of the procedure.

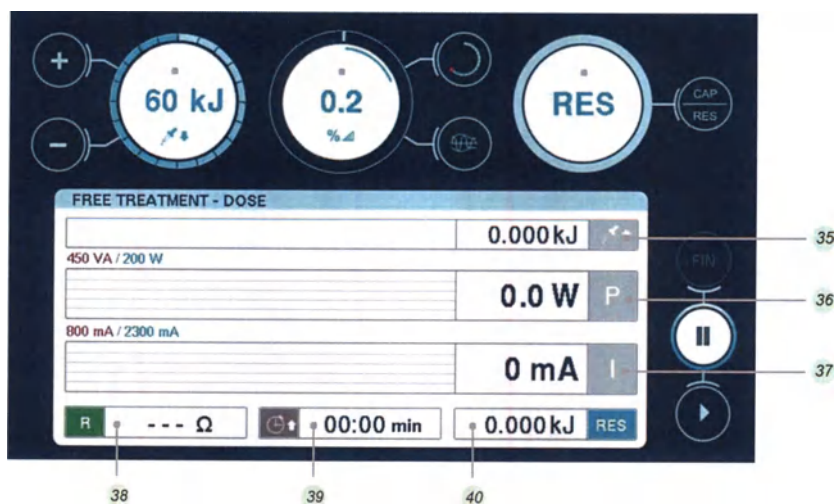
USE OF PROCEDURES

To start or end a procedure:

- Press the "START" button (31) on the device, or the "START / STOP" button on the control panel (14), or set the encoder (5) to the "Ready for contact" state.
- When an electrode contact with a patient is detected, the icon flashes. The countdown of the procedure has not yet begun.
- When the device detects that the electrode is touching the patient, it switches to the status "Procedure in progress", displays the selected parameters of the procedure and starts the countdown of the procedure. The icon indicating that the device has detected contact with the patient's body disappears.
- If the contact between the electrode and the patient disappears during the procedure, the device detects this and returns to the "Ready for contact" state. At the same time, the countdown of the procedure is suspended. When the electrode comes into contact with the patient again, the device switches to the "Procedure is in progress" state.
- To stop the procedure, press the PAUSE button (30) on the device, or the START / STOP button (14) on the control panel, or use the encoder (5).
- To end the procedure and reset the total time of the procedure and the energy transferred to the current point in time, press the END button (29). When you click on this button, a pop-up window appears with a brief description of the various steps in this procedure.
- Press the HOME button (15) to return to the main screen of the unit.

SAVED PROCEDURES

Press the button (23) to access the stored procedures.

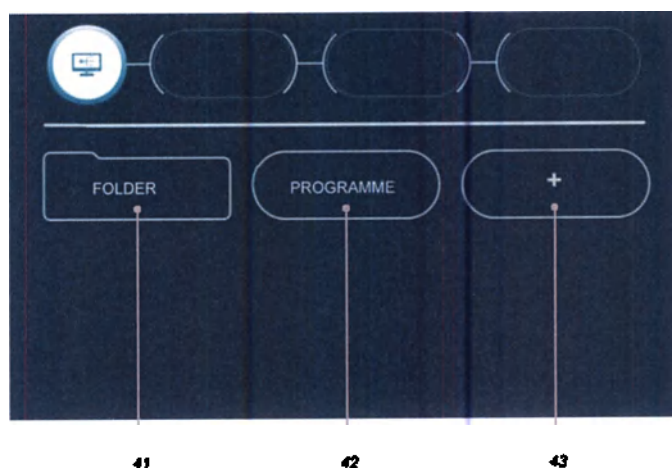- 41 Folder selection button.
- 42 Program selection button.
- 43 Button to add a folder or program.

Above the dividing line on the screen are navigation buttons. They tell the user which folder is currently

open. There are three possible levels. In the first two, both folders and programs are allowed, and the third level is intended only for programs.


The button highlighted in the navigation area indicates to the user at what stage of the selection of the stored procedure he is.

To change the type of procedure, the user must use the buttons above the dividing line on the display.

The buttons below the dividing line are buttons for selecting or entering stored procedures.

Creating Procedures

To create a procedure, press the button (43), after which the following screen will appear:

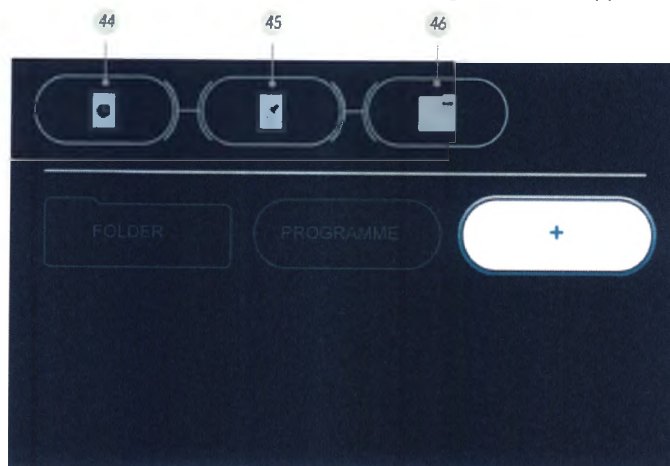44. Creating a time-based program.

45. Creating a program based on the amount of energy transferred.

46. Creating a folder.

After pressing any of these buttons, enter the name of the new item and click on the confirmation button to continue.

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
[Handwritten signature]

- * To create programs based on the countdown, the device uses the following screen (44):

- 47. The choice of the mode of the stage of the procedure (in this case, CAP).
 - 48. The choice of the mode of the stage of the procedure (in this case, RES).
 - 49. . Timing of the stage of the procedure. The minimum time is 1 minute and the maximum is 99 minutes. Time can be increased or decreased at intervals of one minute.
 - 50. The choice of the output signal (in this case, the signal is modulated).
 - 51. . The choice of the output signal (in this case, the signal is not modulated).
 - 52. Selecting an IAS value. IAS values can be increased or decreased at intervals of 1.
 - 53. . The choice of the limiting value of the output power. The minimum value is 5%, and the maximum is 90%. The value from 5% to 30% can be increased or decreased with an interval of 5%. The value from 30% to 90% can be increased or decreased with an interval of 10%.
 - 54. Select the previous step of the procedure.
 - 55. . Select the next step in the procedure.
 - 56. . Deletion of a procedure step.
 - 57. Increase the selected value.
 - 58. Decrease the selected value.
 - 59. Confirmation of the program.
- After pressing the program confirmation button (59), the unit will return to the selection screen and a new program will appear.

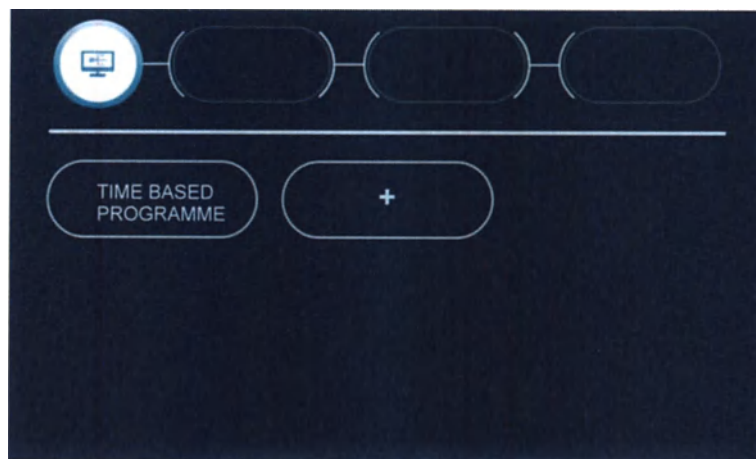
* IAS = INDIBA analogical scale


- To create programs based on the amount of energy transferred, the device uses the following screen (45):


47. The choice of stage mode (in this case, CAP).
48. . The choice of stage mode (in this case, RES).
49. . The choice of the amount of energy transferred to the stage of the procedure (the minimum value is 1 kJ, and the maximum is 999 kJ, the change step is 1 kJ). The energy value can be increased or decreased by 10 kJ by pressing and holding the button for determining the amount of energy transferred.
50. . The choice of the output signal (in this case, the signal is modulated)
51. . The choice of the output signal (in this case, the signal is not modulated).
52. Selecting an IAS value. The IAS value can be increased or decreased with an interval of 1 point.
53. . The choice of the limiting value of the output power. The minimum value is 5%, and the maximum is 90%. The value from 5% to 30% can be increased or decreased with an interval of 5%. The value from 30% to 90% can be increased or decreased with an interval of 10%.
54. . Select the previous step of the procedure.
55. . Select the next step in the procedure.
56. . Deletion of a procedure step.
57. . Increase the selected value.
58. Decrease the selected value.
59. Confirmation of the program. After pressing the program confirmation button (59), the unit will return to the selection screen and a new program will appear.

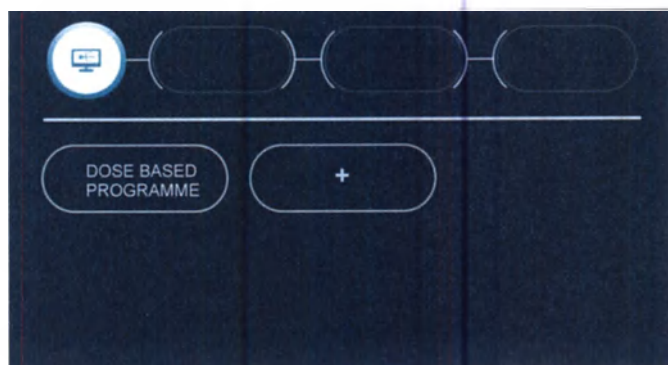When creating a folder (46), the device displays your current location and adds the created folder:

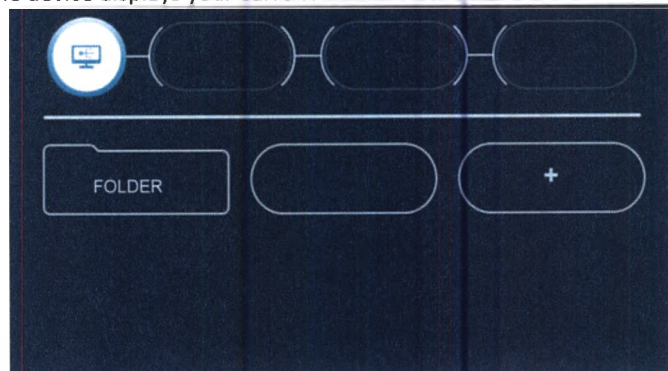
Click on the created folder to go to a new level in the structure, where you can create a new program or folder:

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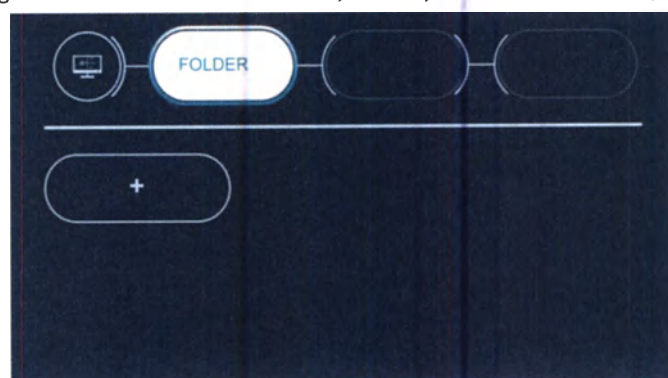Click on the created folder to go to a new level in the structure, where you can create a new program or folder:

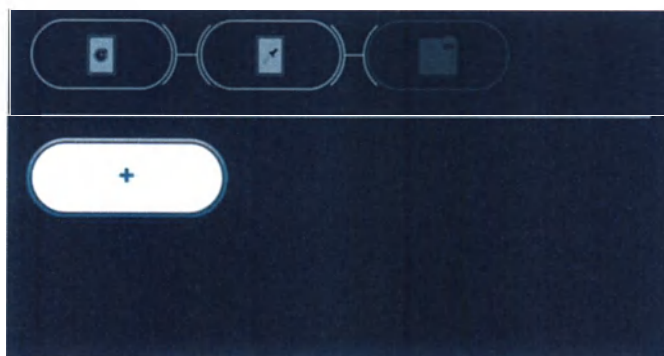Folder management

To manage the folder, hold down the corresponding button for 5 seconds, after which the following options will appear:

60

61

60. . Edit the name of the folder.

61. . Delete the folder (this is possible only if the folder is empty).

Note: The predefined INDIBA folders cannot be changed, deleted, moved or duplicated.

Program management

To control the program, hold down the corresponding button for 5 seconds, after which the following options will appear:

60

61

62

63



- 60. . Change the program.
- 61. . Uninstall the program.
- 62. . Move the program to a folder or subfolder.
- 63. . Copy the program to a folder or subfolder.

Note: The predefined INDIBA folders cannot be changed, deleted, moved or duplicated. They can only be duplicated for use as a template for a new saved program.



Change folders and programs

To exit the program without saving the changes, press the button with the arrow (16). If you save the changes, a message appears on the screen stating that the changes were successfully made:

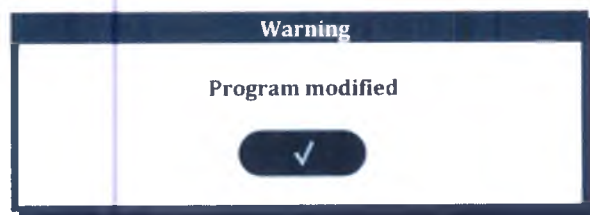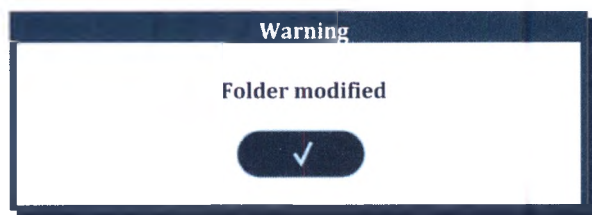
Delete folders and programs

To exit the program without saving the changes, press the button with the arrow (16). If you save the changes, a message appears on the screen stating that the changes were successfully made:

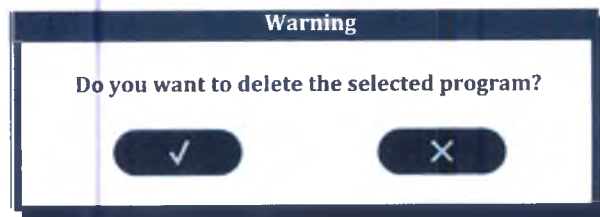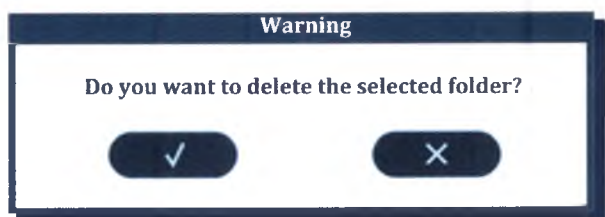If the folder cannot be deleted, the following message appears on the screen stating that the folder cannot be deleted because it contains data:

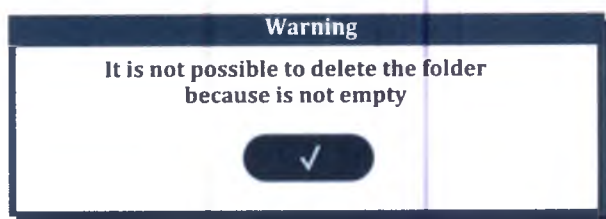
Note: Before deleting, make sure the folder is empty.

Moving Programs

Open folders and subfolders until you are where you want to move the program (the text of the item that will be moved will turn green). Click on an item to confirm that you want to move the program.


A message will appear on the screen informing the user that the program has been correctly moved:

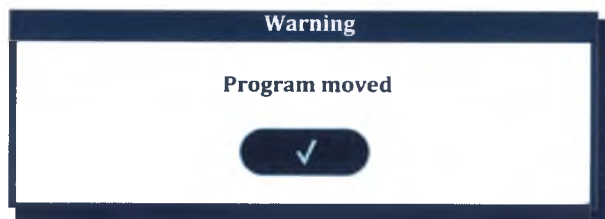Copying Programs

Open folders and subfolders until you are where you want to move the program (the text of the item that will be moved will turn green). Click on an item to confirm that you want to copy the program.


A message will appear on the screen informing the user that the program has been correctly copied:

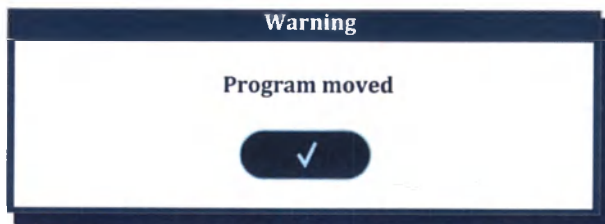Loading saved procedures

Using the buttons under the dividing line of the screen, open folders and subfolders until you find the appropriate program.

When you click the buttons, they appear above the dividing line.

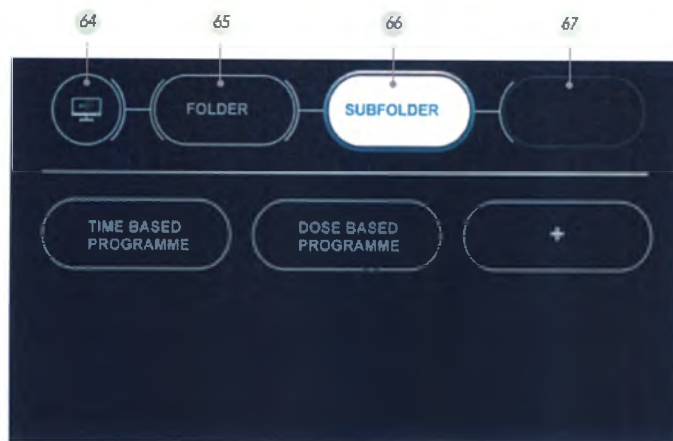64. Screen saved procedures.


65. . The selected folder.

66. . Selected subfolder. Current position.

67. The folder or program to be selected.

If the user has not found the necessary program, he can go back using the buttons 64, 65 and 66, or use the navigation button to return to the previous screen (16).

After selecting a program, the device shows the various stages of the saved program that will be used during the procedure:

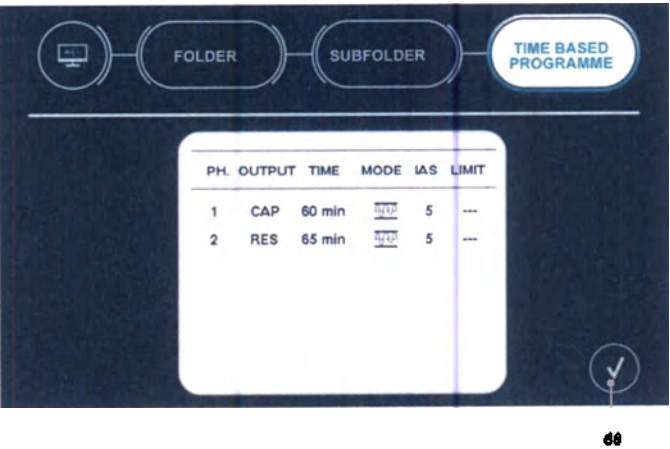68. Downloading the program.

If the user made a mistake when choosing a program, he can return to the previous screen only using the navigation button (16).
When the button (68) is pressed, the program displays the screen of the selected procedure.

INDIBA PREPARED PROCEDURES

In the stored procedures (23) you can find the predefined procedures created by INDIBA. Such folders and / or programs are identified by the following icon:

Settings menu

Press the "Settings" button (17) to go to the settings menu.

You can use the buttons on this screen to:

- 69. Change the volume level of the device.
- 70. Change the sound emitted by the device when the button is pressed.
- 71. . Select a language.
- 72. . Turn off the control of the output power value or set the power limit value: from 5% to 30% (5% interval) and from 30% to 90% (10% interval).
- 73. Select a value (Ω or mA) to show its change during the procedure on the output current graph.
- 74. Change the channel number used for wireless communication with the control panel.
- 75. Change the heating time of the return flat electrode.
- 76. . Turn on or off the sound of contact detection.
- 77. . Decrease the value of the selected parameter.
- 78. Increase the value of the selected parameter.
- 79. Save the changes.

Monitoring program

The button of the monitoring program (18) enables or disables the wireless connection to the remote computer to transmit the values of the parameters measured by the device during the procedure.

Additional information on how the monitoring program works is given in Appendix 1.

MONITORING PROGRAM BUTTON INDICATOR COLOR DESIGNATION

The color of the button indicates the communication status. When the monitoring program icon:

- It is white and does not blink; the Wi-Fi connection is disconnected.
- White and blinking; Wi-Fi connection is on and waiting for data connection.
- It is blue and flashing; a data connection has been established, but data is not sent to the monitoring program.
- It is blue and does not blink, the connection is completed, and the device sends data for storage.

Flat return electrode heating

Note: This function is only available if the unit is mounted on a trolley rack. Before activating this function, the user must connect the device to the trolley rack using the connecting cable and make sure that the selected voltage value on the trolley rack matches the existing mains voltage.

The return plate heating button (19) is used to enable or disable the return plate heating.

The heating time can be changed in the "Settings" menu (17).

The color of the button indicates the status:


- White and does not blink - heating is off.
- Yellow and does not blink; heating is off.
- Yellow and flashing; heating will soon be turned off.

If you press the button during heating, the countdown will be interrupted.

Information on using the device

When you press the information button on the use of the device (20), a message appears on the screen with the following brief data:

- The number of procedures performed on the device.
- Total time of procedures performed on the device.
- Total time of procedures performed on the device in CAP mode.
- Total time of procedures performed on the unit in RES mode.
- The amount of time remaining until the next maintenance.

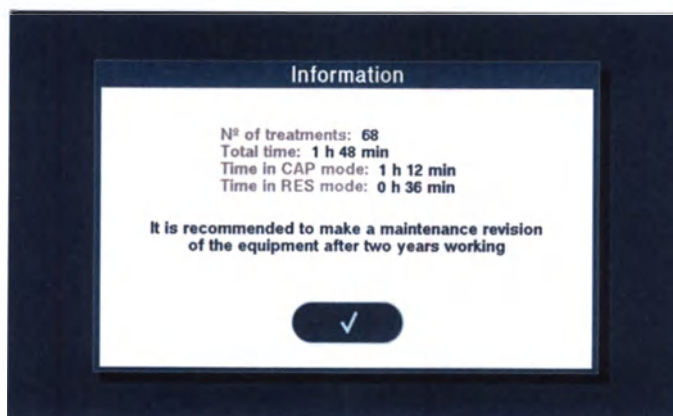DEVICE USE BUTTON

The information button on the use of the device (20) indicates to the user when it is necessary to perform maintenance on the device.

- The button is white and does not blink. This means that there is still time until the next service date for the machine.
- The button is blue and does not blink. This means that the device must be serviced in the near future.
- The button is red and does not blink. This means that the service day of the device has already passed.

The user must contact their dealer or INDIBA customer support to service the machine.

Icons

During the procedure, the following yellow icons may flash at the bottom of the screen:


80. The selected output value cannot be reached because the resistance value detected by the device is too large.

81. The detected capacitive electrode is too large for this unit.

82. . This icon appears in capacitive mode (CAP) when the device adjusts its internal parameters to adapt to the detected electrode and optimize performance.

83. The device detected contact of the electrode with the patient.

If the device cannot detect the connection of the electrode handles, it will stop the procedure and display a pop-up message on the screen for the user:

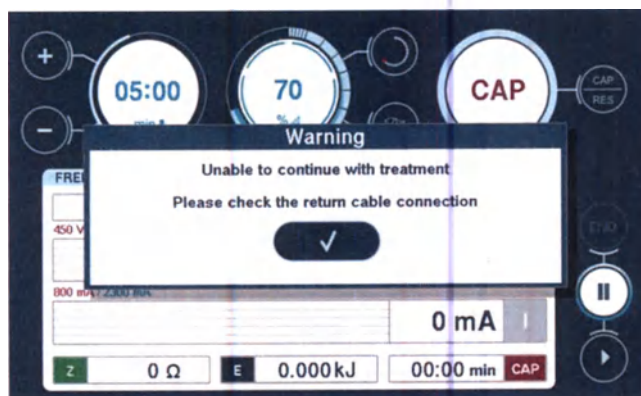The user needs to check the cable connection indicated in the message.

REMOTE CONTROL

The control panel comes with one tablet battery. If the unit does not respond to the control panel, it may be necessary to replace the battery. Remove the battery from the control panel if it will not be used for a long time.

When the battery in the control panel runs out, a solid yellow icon appears on the display of the device.

84. The battery in the control panel is empty, it must be replaced as soon as possible. Follow these instructions to replace the battery:

- Remove the battery cover.
- Remove the discharged battery.
- Insert a new battery into the battery compartment.
- Replace the battery cover.
- Dispose of the old battery at a recycling point.

After the device has reached the end of its useful life, remove the battery from the control panel and return it to a recycling center for environmentally friendly disposal before disposing of it.

Product Elements

ELECTRODE TYPES

The following types of electrodes are used in the apparatus:

- Returnable
- Capacitive
- Resistive

RETURN ELECTRODES

Two types of electrodes act as a return electrode:


- Return cylindrical electrode.
- Return flat electrode.

A return flat electrode should be used when the diameter of the active electrode (capacitive or resistive) is greater than 40 mm. When using an active electrode (capacitive or resistive) with a smaller diameter, a flat or cylindrical return electrode can be used, depending on the patient's sensations, the area to be treated, etc.

Notes: It is necessary to use a return electrode (flat or cylindrical) so that the current can circulate.


The return flat electrode is flexible, allowing you to better adapt to the anatomy of the patient's body. However, the following bending limits of the return flat electrode must be considered:

- In the transverse direction (A), the return flat electrode cannot be bent more than 30° horizontally.
- In the longitudinal direction, the return flat electrode cannot be bent more than 90° provided that pressure is applied to the extreme parts of the return flat electrode (B).
- NEVER bend the return flat electrode in the longitudinal direction, applying pressure at points close to the bend line (C), as this will deform the return flat electrode for a long time, making it impossible to return to its original shape.

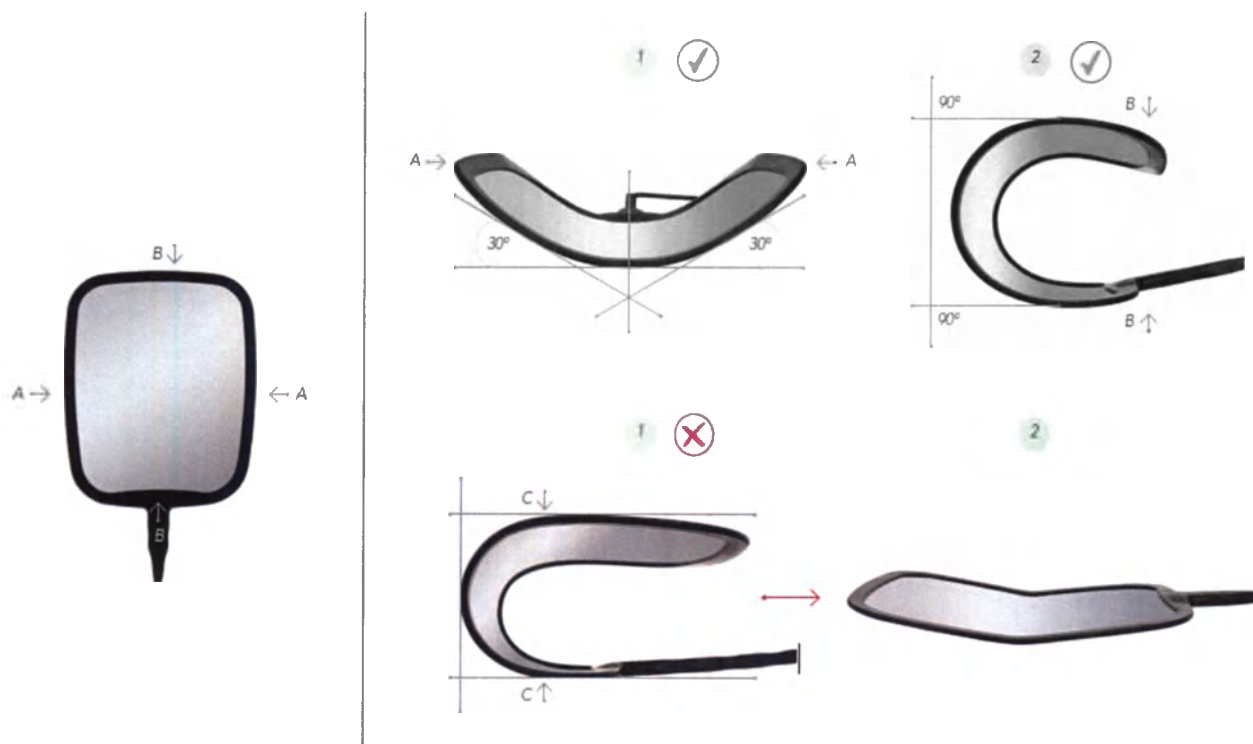HANDLE FOR ELECTRODES CAP


Only capacitive electrodes with a shaft diameter of 3 mm and a groove on the shaft can be used in the product. This groove helps to fix the electrode in the required position on the electrode handle.

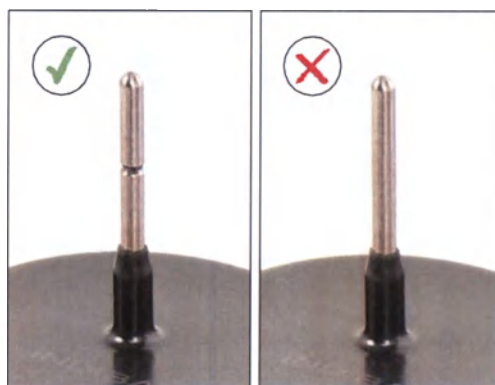In Fig. Figures 1-4 show the steps of fixing a capacitive electrode on a handle for electrodes:


85. Insert the protective sleeve into the electrode handle.
86. Insert the electrode shaft into the hole until it stops.

1. Release the protective sleeve to return it to its original position.
2. Gently pull the electrode until you hear a click. This indicates that the groove on the shaft is correctly installed in the handle for the electrodes.

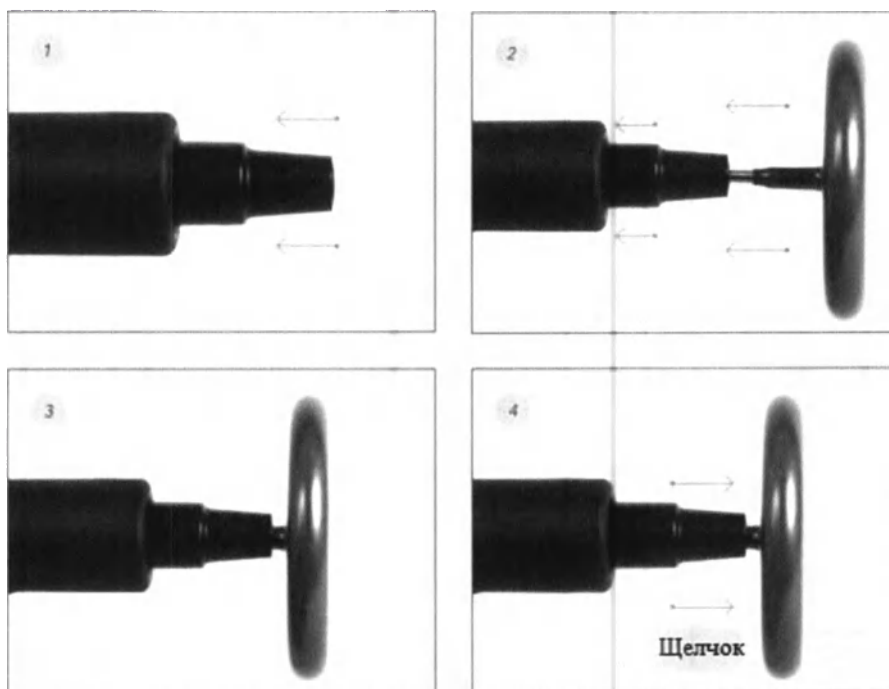Figures 1-3 show the steps for removing a capacitive electrode from the electrode handle:

1. Insert the protective sleeve into the electrode handle.
2. Pull the capacitive electrode out of the electrode handle.
3. Release the protective sleeve to return it to its original position..

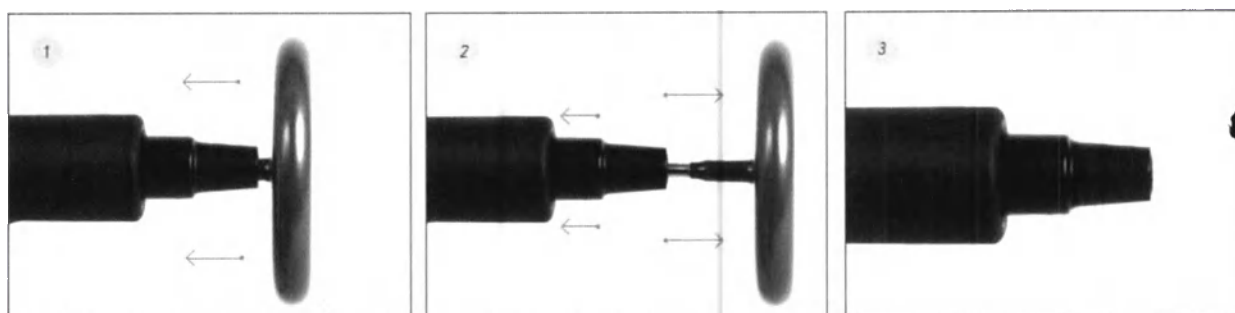IMPORTANT: To insert, pull out or replace an electrode, the procedure must be paused. Do not start the procedure unless you hear a click indicating that the electrode is locked. Otherwise, the electrode may fall during the procedure and break.


HANDLE FOR RES ELECTRODES

Only resistive electrodes with a shaft with a diameter of 6 mm and a groove on the shaft can be used in the apparatus. This groove helps to fix the electrode in the required position on the electrode handle.


Figures 1-4 show the steps of fixing a resistive electrode on a handle for electrodes:

1. 1. Insert the protective sleeve into the electrode handle
2. 2. Insert the electrode shaft into the hole until it stops.
3. 3. Release the protective sleeve to return it to its original position.
4. 4. Gently pull the electrode until you hear a click. This indicates that the groove on the shaft is correctly installed in the handle for the electrodes.

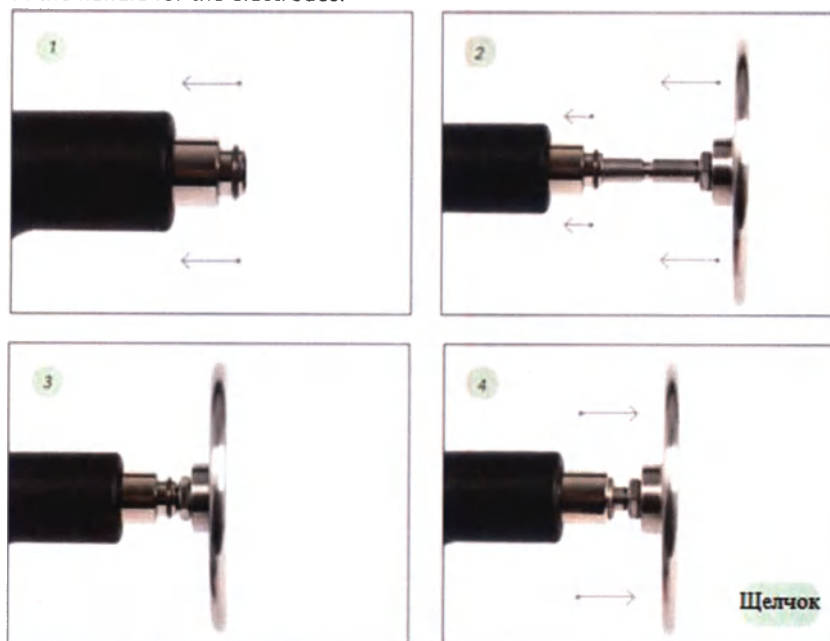In Fig. 1-3 shows the steps for removing a capacitive electrode from the handle for the electrode:

1. 1. Insert the protective sleeve into the electrode handle.
2. 2. Pull the resistive electrode out of the electrode handle.
3. 3. Release the protective sleeve to return it to its original position.

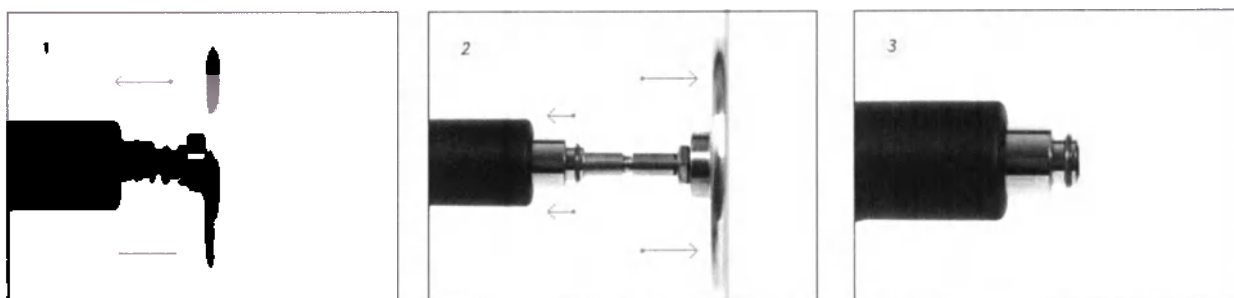ВАЖНО: To insert, pull out or replace the electrode, the procedure must be put on PAUSE. Do not start the procedure unless you hear a click indicating that the electrode is locked. Otherwise, the electrode may fall during the procedure and break.

MANUAL


To start the procedure:

- Correctly fix the return electrode on the patient's body.
- Correctly install the active electrode on the electrode handle.
- Apply a conductive cream to the area of the patient's body to be treated and to the return electrode.



PLACING AN ACTIVE ELECTRODE AND RETURN ELECTRODE ON THE BODY OF THE PATIENT

The figure below shows the recommended placement of the active electrode and the return electrode based on which area of the body is to be processed.


1. Skin remodeling

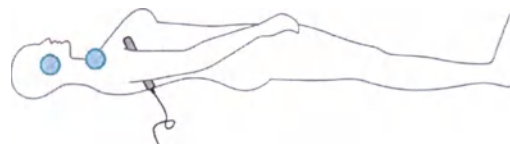2. Body contouring

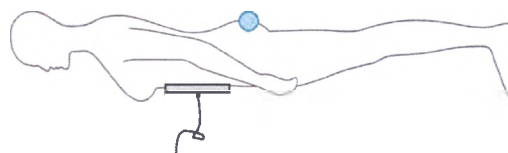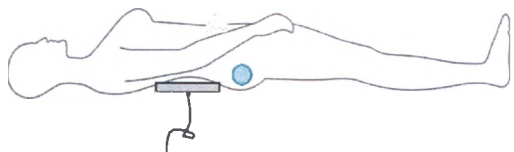
3. Anti-cellulite programs

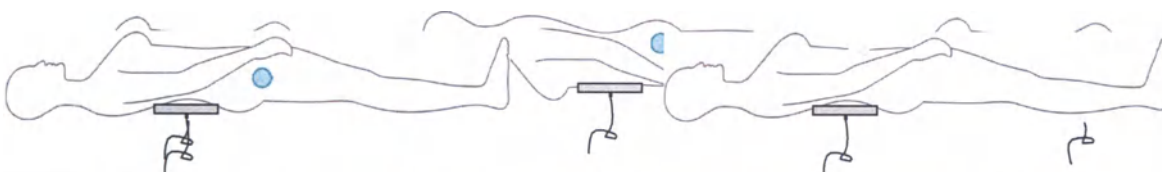
4. Postoperative recovery

Liposuction

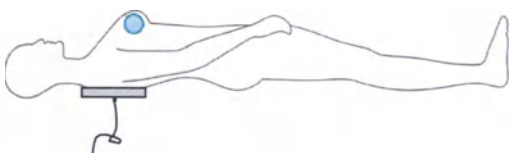Abdominoplasty

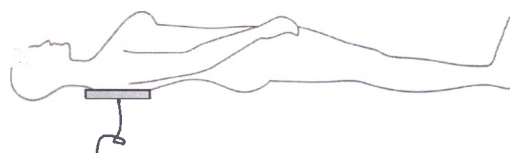Mammoplasty

Face lift


IMPORTANT: Never use the unit if the return electrode has not been properly fixed. Failure to do so may result in burns to the patient.

IMPORTANT: Do not use electrodes on the patient's body until you have fully studied this instruction manual and completed the training course provided by the manufacturer when purchasing the device.

GUIDELINES FOR USING ACTIVE ELECTRODES


Capacitive Electrodes (CAP)

The electrode should always come in contact with the skin. During the procedure, the electrode should move. Use circular, linear, or both movements in accordance with recommendations or treatment protocols.

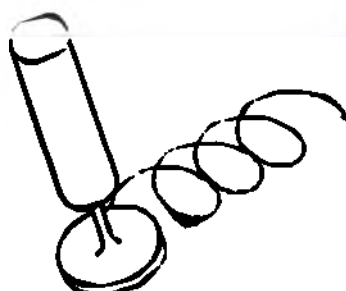Apply moderate pressure and move forward in a wide circular motion.

The electrode should fit snugly against the skin. Using the electrode at an angle at which only partial contact with the skin is observed can cause premature wear of the edges of the electrode, increase the risk of burns and reduce the effectiveness of the procedure.

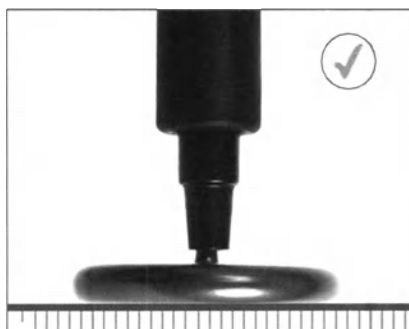Capacitive electrode (CAP).

Д To obtain optimal treatment results, three key points must be observed:

Pressure

The applied pressure should be more or less the same as if using a large eraser. Under these conditions, the minimum electrical resistance is observed and the best physical contact is ensured.

Speed

In each circle, the electrode should move a distance equal to 50% of its diameter. The electrode should move continuously at a speed of approximately one rotation per second.

Temperature

Choose a comfortable procedure temperature. The patient should feel warm, but this sensation should not cause discomfort.

The ratio of the diameter of the capacitive electrode (CAP) to the minimum and maximum effective time of the procedure for an area of 100 cm².

Diameter (mm)	Min time (min)	Max. time (min)
30	10	20
40	5	10
55	2:30	5
65	1:30	3
80	1	2

The duration shown is approximate. The exact value is selected.

This table shows the relationship between the diameter of capacitive electrodes (CAP) and the procedure time.

IMPORTANT: If damage to the capacitive electrode is detected, it must be replaced. Failure to do so may result in minor burns to the patient's skin.

Resistive Electrodes (RES)

The electrode should always come in contact with the skin. During the procedure, the electrode should move. Use circular, linear, or both movements in accordance with recommendations or treatment protocols.

Apply moderate pressure and move forward in a wide circular motion..

The electrode should fit snugly against the skin. Using the electrode at an angle at which only partial contact with the skin is observed can cause premature wear of the edges of the electrode, increase the risk of burns and reduce the effectiveness of the procedure.

Resistive electrode (RES).

To obtain optimal treatment results, three key points must be observed:

Pressure

The applied pressure should be more or less the same as if using a large eraser. Under these conditions, the minimum electrical resistance is observed and the best physical contact is ensured.

Speed

In each circle, the electrode should move a distance equal to 50% of its diameter. The electrode should move continuously at a speed of approximately one rotation per second.

Temperature

Choose a comfortable procedure temperature. The patient should feel warm, but this sensation should not cause discomfort.

NOTES

- Before starting the procedure in resistive mode, it is recommended that you first perform the procedure in capacitive mode.
- Procedure time: set at the discretion of the specialist, based on the parameters of the procedure and processing area. The duration of the procedure can vary from 20 minutes for procedures performed on the face, and up to 50 minutes for procedures performed on the body.
- The effectiveness of the procedure is determined by two main factors: the duration and level of output power.
- Hard hair must be removed from the treatment area.
- Resistive mode can sometimes cause internal discomfort to patients in the treatment area. Reduce the output to a level where the patient feels warm but not uncomfortable.

CABLE CONNECTIONS USER GUIDE

1. Cables should be stored separately from each other, take a natural path created by the weight of the cable itself.
2. Do not overtighten or bend the cables. Cables should not be entangled or twisted, as this can lead to premature damage.

3. Disconnect the cable from the device as shown in fig.

4. Do not pull the cable to disconnect it.

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97
2020

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[Signature]

Electroconductive creams manufactured by INDIBA

Electroconductive creams manufactured by INDIBA were specially developed for devices manufactured by INDIBA.

Apply an electrically conductive cream to the skin of the treated area (layer thickness should be at least 1 mm). Add an electrically conductive cream during the procedure if you notice that it has lost its original properties.

Apply a conductive cream to the return electrode before starting the procedure.

INDIBA electrically conductive creams have two main benefits for the patient:

- They reduce the electrical resistance of the electrode-skin contact, thereby optimizing the flow of current through the body.
- They minimize friction between the skin and the electrode, allowing the electrode to glide smoothly over the skin. The package indicates the composition of the cream, batch and expiration date.

Product cleaning and disinfection

CLEANING THE UNIT

- Always keep the unit in good condition. Periodically clean it from dust and other impurities with a soft, slightly damp cloth. Do not allow water droplets to fall into the connectors of the unit. Do not use aggressive or oxidizing substances such as bleach or ammonia.

Cleaning and disinfection of active electrodes, return electrodes and electrode handles

Before and after each procedure

- Always clean active, return electrodes and electrode handles before and after each use.
- Use a lint-free cloth or towel to clean excess conductive cream from active electrodes, return electrodes, and electrode handles.
- Then, active, return electrodes and electrode handles must be disinfected with a disinfectant wipe suitable for non-invasive medical instruments, or a mild soap solution and water. Dry before use.
- It is recommended to use a new cloth to clean each item.

Capacitive (CAP) electrodes are not intended for sterilization. Sterilization can damage the coating.

After each use, the massage nozzle should be washed with soapy water and then dried..

When carrying out gynecological procedures:

Longitudinal electrode RES:

For individual use only. Clean the electrode twice with soapy water. Then apply a disinfectant that is intended for professional use (phenol or glutaraldehyde 2%) for 30 minutes, and then rinse with water.

Longitudinal CAP electrode:

For individual use only. Clean the electrode twice with soapy water. Disinfect the electrode handles and return electrodes using the equipment used in the hospital. If in doubt, contact the manufacturer.

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Repair and maintenance

Daily inspection

Before using the device, make sure that the return and active electrodes, electrode handles are in good condition: the case is intact, the cable coating is not damaged. Pay particular attention to the coating on the capacitive electrodes; it must not be damaged.

Do not use damaged or defective product elements. If repairs or replacements are necessary, contact the manufacturer or your authorized dealer.

Product repairs may only be carried out by authorized specialists or service centers. Contact the manufacturer for information on the nearest authorized service center.

Periodic inspections by a service engineer should be carried out at least once every two years.

IMPORTANT: If any defect or malfunction is detected, DO NOT use this product, but notify a qualified technician. To avoid electric shock, DO NOT disassemble this product by yourself.

Marking and packaging

Marking

Marking of product elements may contain the following applicable parameters, symbols and designations:

Product name;
 Name of product item;
 Input electrical parameters: voltage, frequency of power supply, power;
 Class of protection against electric shock;
 Parameters of electromagnetic compatibility; Information about the authorized representative of the manufacturer in the Russian Federation;
 Parameters of electromagnetic compatibility; Information about the authorized representative of the manufacturer in the Russian Federation;
 Fuse information;
 Output data: power output in CAP and RES mode, frequency

Symbol "Serial number" with the serial number;

Symbol "Catalog number" indicating the catalog number;

The "Use To" symbol. Indicates the date after which the product should not be used;

Symbol "Part Number". Indicates the medical device number according to the manufacturer's catalog;

Symbol "Date of manufacture" indicating the date of manufacture;

Symbol "Information about the manufacturer" with the name of the manufacturer and address (optional: telephone and fax);

Symbol "Working parts of type BF without protection against defibrillator discharge";

Symbol "WEEE directive - Do not dispose of as normal household waste";

Symbol "CE mark of conformity" with (or without) indicating the number of the notifying authority;

Without the use of parabens;

Refer to the product manual;

After use, the packaging should be thrown into the bin;

Product packaging is made of plastic and is amenable to industrial processing;

Shelf life after opening the package;

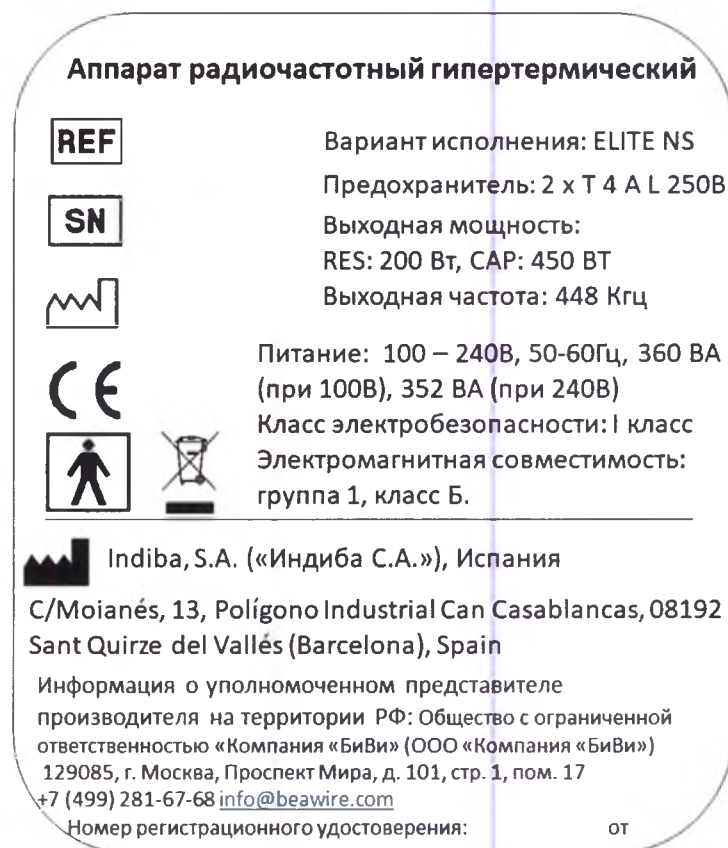
Fig. 45 Design of marking the main ELITE NS unit (for the main RGN unit - similarly)

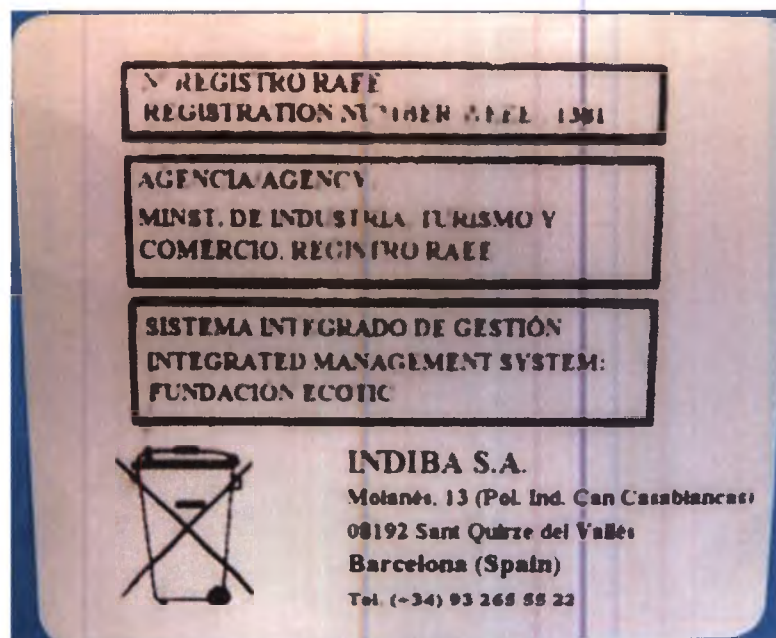
Fig. 46 Additional marking of the device, not applicable for Russia

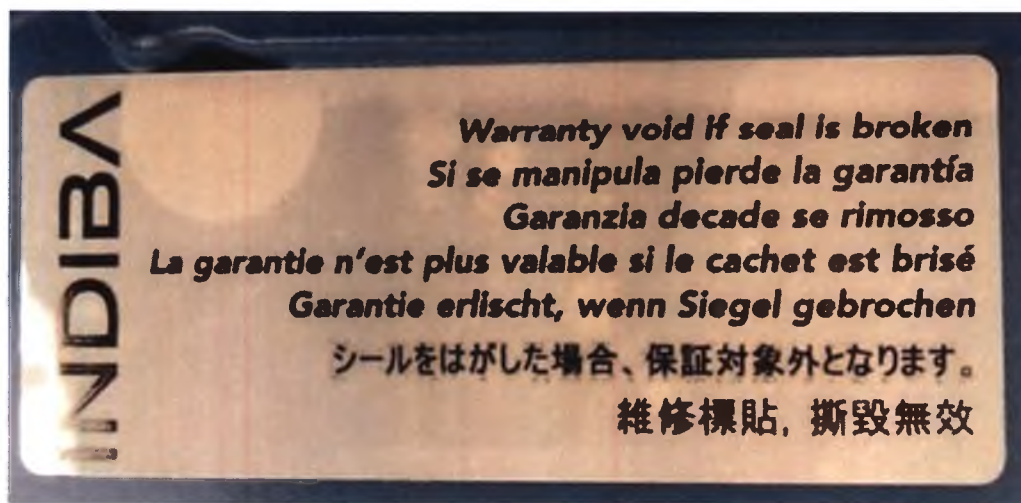
Fig. 47 Warranty seal of the device "Warranty is void if the seal is damaged"

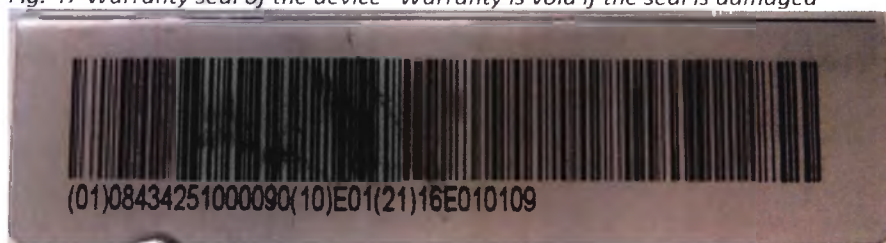
Fig. 48 Marking on the back of the machine containing a unique product identification number

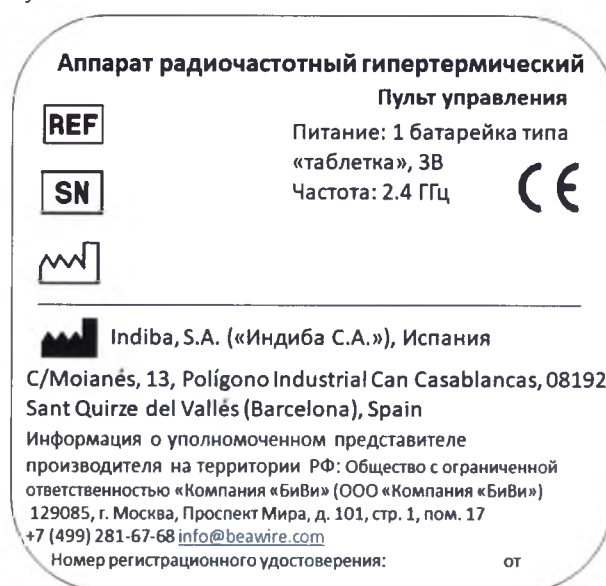
Fig. 49 Remote control labeling project

Fig. 50 Marking in the packaging of the apparatus on the procedure for removing the product from packaging

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Fig. 52 Additional labeling of Anticellulite Cream for Russian user

INGREDIENTS
INGREDIENTES
COMPOSITION
INGREDIENTI
INGREDIENTES
INHALTSSTOFFE
INGREDIENTEN

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Aqua
PEG-6 Stearate
Paraffinum liquidum
Hydrogenated Palm/palm kernel oil PEG-6 esters
Caprylic/capric Triglyceride
Glycol Stearate
PEG-32 Stearate
Sodium Sulphate
PEG-75 Lanolin
Benzyl Alcohol
Triethanolamine
Salicylic Acid
Glycenn
Sorbic Acid

CA INSTRUCTIONS FOR USE: Apply to the area to be treated by spreading a sufficient amount of cream so that the average thickness deposited on the skin is not less than 1 mm. Add more cream during the treatment when it is noticed that it no longer has its initial smoothness and its density increases, making hand movements difficult. Remove the remaining cream after treatment.

WARNING: Only Proionic Creams have been designed and tested for use with INDIBA Proionic Systems.

Use of other gels or creams invalidates the warranty of the equipment and may cause harm to the person under treatment.

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Copyright © 2015 Indiba S.A. All rights reserved. Made in Spain

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Fig. 54 Additional labeling of conductive face cream for Russian user

INDIBA
PROIONIC CREAMS

INDIBA
PROIONIC CREAMS

CONDUCTIVE CREAM FOR
CAPACITIVE & RESISTIVE TREATMENT

1000 ml. e 35 fl oz

INGREDIENTS
INGREDIENTES
COMPOSITION
INGREDIENTI
INGREDIENTES
INHALTSSTOFFE
INGREDIËNTEN

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Aqua
PEG-6 Stearate
Paraffinum liquidum
Hydrogenated Palm/palm kernel oil PEG-6 esters
Caprylic/capric Triglyceride
Glycol Stearate
PEG-32 Stearate
Sodium Sulphate
PEG-75 Lanolin
Parfûm
Benzyl Alcohol
Triethanolamine
Salicylic Acid
Glycerin
Sorbic Acid

Copyright 2015 Indiba S.A. All rights reserved. Made in Spain

GB INSTRUCTIONS FOR USE: Apply to the area to be treated by spreading a sufficient amount of cream so that the average thickness deposited on the skin is not less than 1 mm. Add more cream during the treatment when it is noticed that it no longer has its initial smoothness and its density increases, making hand movements difficult. Remove the remaining cream after treatment.

WARNING: Only Proionic creams have been designed and tested for use with INDIBA Proionic Systems.

Use of other gels or creams invalidates the warranty of the equipment and may cause harm to the patient.

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Fig. 55 Original marking of conductive Care Cream

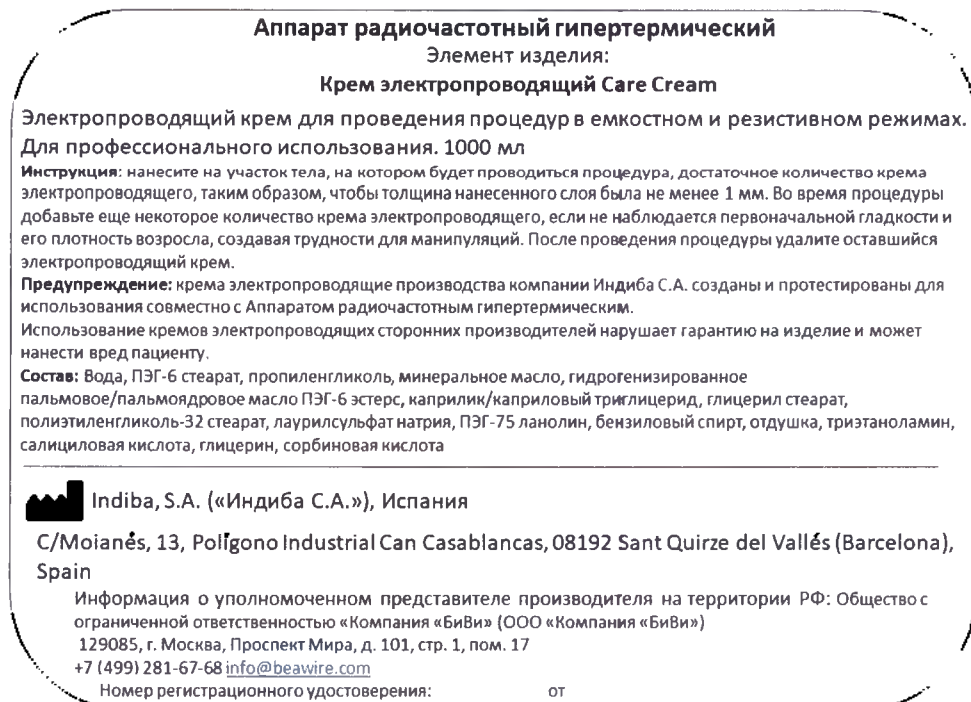
Fig. 56 Additional marking of cream of electrically conductive Care Cream for a Russian-speaking user

Packaging

The main unit is delivered in a cardboard box measuring 590mm * 610mm * 380mm, ± 10 mm. When it is opened, the user will see an explanatory marking with information on the correct removal of the device from the shipping container. Also, a sheet of foamed polyethylene is placed on top as an additional protection.

Fig. 57 Transport packaging of the main unit

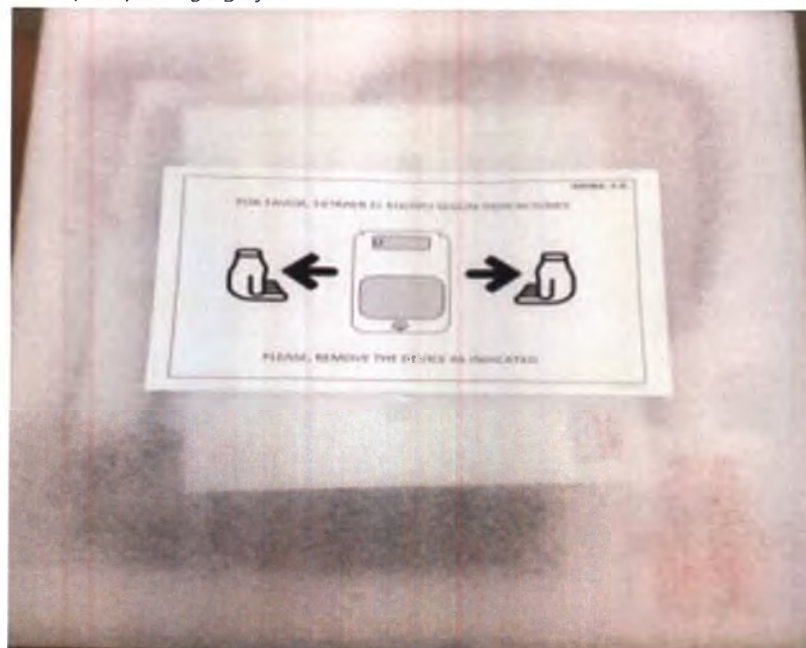
Fig. 58 Explanatory marking on the correct removal of the device from the packaging

When removing the protective sheet of foamed polyethylene, the user will see elements of the composition of the product, each of which is in its cell form of foamed polyethylene.

Fig. 59 Product elements (active electrodes (RES and CAP) are under the return flat electrode)

After extracting the form with the elements of the product, the user will find the main unit in a non-woven bag, which is equipped with shock protection from foamed polyethylene.

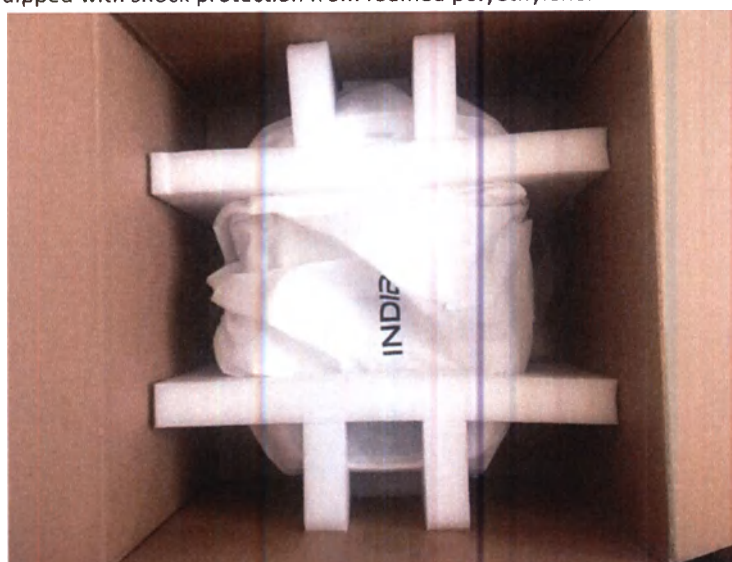
Fig. 60 Packing the main unit

Fig. 61 Marking of the transport packaging of the device using an example of an embodiment Radio-frequency hyperthermic device ELITE NS (for another embodiment - similarly). The marking indicates that the product is fragile (FRAGILE), as well as that made in Spain (made in Spain).

The trolley is delivered in a cardboard box measuring 600mm * 400mm * 900mm, ± 10 mm, tied with packing tape.

Fig. 62 Transport packaging rack trolley

When opening the transport packaging, the user will find a trolley-rack in a plastic case, on the edges of which a foam polyethylene is worn to protect it from shock. The trolley stands on a wooden frame.

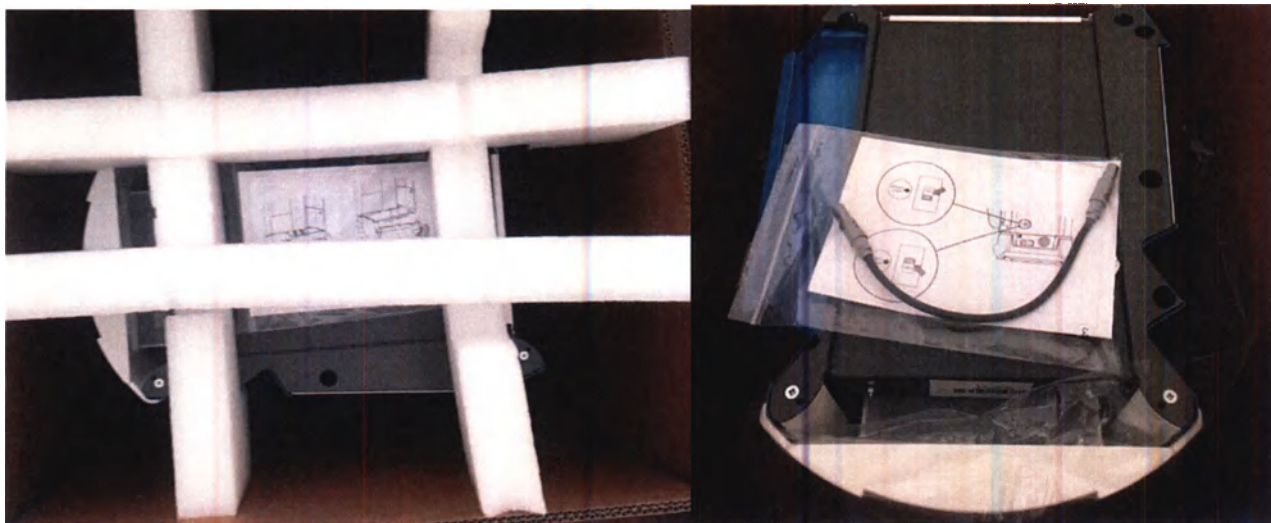
Fig. 63 Trolley rack with protection from foamed polyethylene, together with the trolley rack graphic material is supplied for turning on and installing the device, a connecting cable and keys for drawers of the trolley rack.

The electrically conductive cream is delivered sealed in a plastic film and a protective cap on the nose:

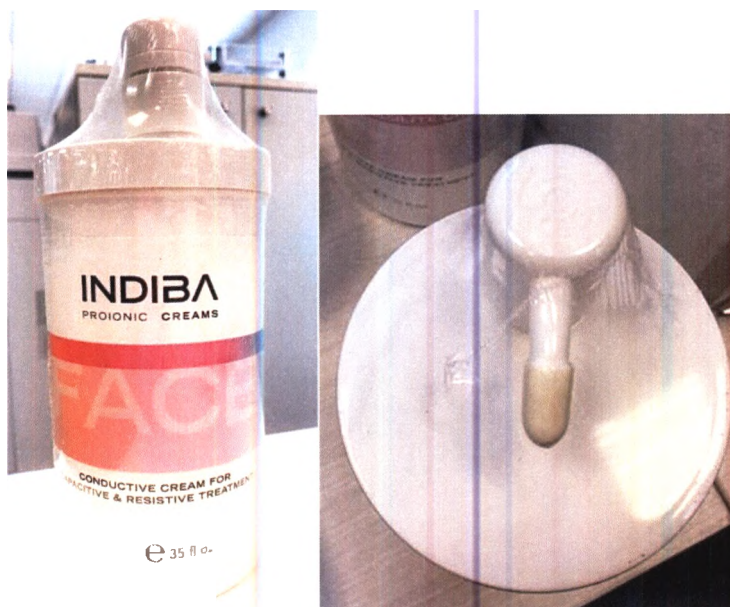
Fig. 64 Conductive Face Cream packed in plastic wrap

Conductive creams can be delivered in a separate transport package of 420mm * 330mm * 250mm $\pm$ 10mm in size or in a transport package of the Welcome Pack:

Fig. 65 Transport packaging cream electrically conductive

Welcome Pack comes in a separate cardboard box measuring 331mm * 365mm * 280mm $\pm$ 10mm. Inside are filled with air-filled plastic bags for protection against shock and damage. This box contains electrically conductive creams (if there is no transport packaging for electrically conductive creams), a bag for electrodes, as well as a separate box with magnetic locks measuring 258mm * 60mm * 295mm $\pm$ 10mm, which contains massage nozzles and a USB drive.

Fig. 66 Welcome Pack Transport Package

Fig. 67 Packing massage tips and a USB stick in a separate box. Graphic material with instructions for installing the massage nozzle is also embedded inside

The PACK GINECO-AESTHETICS kit is supplied in a separate box with magnetic locks (similar to the box described above with the Welcome Pack). This box is placed in a shipping cardboard box measuring 163mm * 423mm * 325mm $\pm$ 10mm.

Fig. 68 Transport packaging for the PACK GINECO-AESTHETICS kit (left) and consumer packaging for the PACK GINECO-AESTHETICS kit (right)

The following designations of the conditions of transportation and storage are available on transport containers:

Top

Fragile, handle with care

Keep away from moisture

Marking of product shipping containers contains:

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Symbol "Information about the manufacturer" (indicating the name of the manufacturer and address).

Technical details

PRODUCT CLASSIFICATION

The product is classified in accordance with international standards and relates to medical electrical equipment with the following characteristics:

- Product electrical safety class: I electrical safety class
- Electromagnetic compatibility: group I, class B
- Type of protection of the working part against leakage currents: working part type BF
- Degree of protection of a cover: IPX0
- Sterilization method: not applicable
- Protection against fire hazards of flammable mixtures of anesthetics with air or with oxygen, or nitrous oxide: no
- Operating mode: long
- Estimated Life: 5 years
- Embedded software: V1.09 from 05/02/2018
- Risk class of firmware: A
- Monitoring program: 0.5R0 of 04/19/2016
- Risk class of firmware: A

COMPLIANCE WITH STANDARDS OF THE RUSSIAN FEDERATION

- GOST ISO 13485-2011 Medical devices. Quality Management Systems. Regulatory system requirements
- GOST ISO 14971-2011 Medical devices. Application of risk management to medical devices
- GOST R ISO 15223-1-2014 Medical devices. Symbols used for labeling on medical devices, labels and accompanying documentation. Part 1. Basic requirements
- GOST R IEC 60601-1-2010 Medical electrical equipment. Part 1. General safety requirements taking into account the basic functional characteristics
- GOST R IEC 60601-1-2-2014 Medical electrical equipment. Part 1-2. General safety requirements, taking into account the basic functional characteristics. Parallel standard. Electromagnetic compatibility. Requirements and Tests
- GOST R 50444-92 Medical devices, apparatus and equipment. General specifications
- GOST R IEC 60601-1-6-2014 Medical electrical equipment. Part 1-6. General safety requirements, taking into account the basic functional characteristics. Additional standard. Serviceability
- GOST R IEC 62304-2013 Medical devices. Software. Life cycle processes
- GOST R ISO / IEC 9126-93 Information technology. Evaluation of software products. Quality characteristics and guidelines for their use
- GOST 28195-89 Software Quality Assessment. General Provisions
- GOST R 51188-98 Information Security. Testing software for computer viruses. Sample manual
- GOST R ISO / IEC 12119-2000 "Information technology. Software packages. Quality Requirements and Testing

- • GOST 31214-2016 Medical devices. Requirements for samples and documentation provided for toxicological, sanitary and chemical tests, sterility and pyrogenicity tests
- • GOST R 52770-2016 Medical devices. Safety requirements. Methods of sanitary-chemical and toxicological tests
- • GOST ISO 10993-1-2011 Medical devices. Assessment of the biological effects of medical devices. Part 1. Assessment and research
- • GOST R ISO 10993-2-2009 Medical devices. Assessment of the biological effects of medical devices. Part 2. Animal Handling Requirements
- • GOST ISO 10993-5-2011 Medical devices. Assessment of the biological effects of medical devices. Part 5. Cytotoxicity assay: invitro methods
- • GOST ISO 10993-10-2011 Medical devices. Assessment of the biological effects of medical devices. Part 10. Studies of irritating and sensitizing effects
- • GOST ISO 10993-12-2015 Medical devices. Assessment of the biological effects of medical devices. Part 12. Sample preparation and control samples
- • GOST 31209-2003 Containers for blood and its components. Chemical and biological safety requirements and test methods
- • MUK 4.1.742-99 Inverse voltammetric measurement of the concentration of zinc, cadmium, lead and copper ions in water
- • GOST 4011-72 Drinking water. Methods for measuring the mass concentration of total iron (with Changes N 1, 2)
- • MUK 4.1.3166-14 Gas chromatographic determination of hexane, heptane, acetaldehyde, acetone, methyl acetate, ethyl acetate, methanol, isopropanol, acrylonitrile, n-propanol, n-propyl acetate, butyl acetate, isobutanol, n-butanol, benzene, toluene, ethyl , m-, o- and p-xylenes, isopropylbenzene, styrene, α -methylstyrene in water and aqueous extracts from materials of various compositions
- • MP 1849-78 Guidelines for the determination of formaldehyde in aqueous extracts and model environments
- • I-880-71 Instruction for the sanitary-chemical research of products made from polymeric and other synthetic materials intended for contact with food products "(ethylene glycol)
- • GOST 31866-2012 Drinking water. Determination of element content by inversion voltammetry
- • GOST 31870-2012 Drinking water. Determination of element content by atomic spectrometry methods (as amended)
- • GN 2.3.3.972-00 Maximum allowable amounts of chemicals released from materials in contact with food
- • MUK 4.1.3169-14 Gas chromatographic determination of dimethyl phthalate, dimethyl terephthalate, diethyl phthalate, dibutyl phthalate, butyl benzyl phthalate, bis (2-ethylhexyl) phthalate and dioctyl phthalate in water and aqueous extracts from materials of various compositions
- • MP 1503-76 Guidelines for the determination of hexamethylene diamine in water during sanitary-chemical studies of polymeric materials used in the food and textile industries
- • MUK 4.1.647-96 Guidelines for the gas chromatographic determination of phenol in water

TRANSPORTATION, STORAGE, OPERATION.

Terms of Use:

- Air temperature: + 10 ° C to + 40 ° C.
(for longitudinal electrode RES and longitudinal electrode CAP, air temperature is from +10 ° C to +42 ° C)
- Relative humidity: 30% - 75%.
- Atmospheric pressure: 70 kPa - 106 kPa.

[Handwritten signature]

Transportation and storage conditions:

- Air temperature: -20 °C - +50 °C.
- Relative humidity: 10% - 100%, non-condensing.
- Atmospheric pressure: 50 kPa - 106 kPa.

Before turning on the device, wait until its temperature becomes equal to room temperature.

WEIGHT-DIMENSIONAL CHARACTERISTICS

Item	Parameter	Value
Block main ELITE NS	Weight, kg, $\pm 1\%$	8.6
	Size, mm, $\pm 1\%$	400x530x150
Main RGN Block	Weight, kg, $\pm 1\%$	8.4
	Size, mm, $\pm 1\%$	400x530x150
Trolley	Weight, kg, $\pm 1\%$	37.72
	Size, mm, $\pm 1\%$	445x325x905
Connecting cable	Weight, kg, $\pm 10\%$	With open drawers, the width increases to 740 mm
	Wire length, mm, $\pm 5\%$	0,025
	Connector Diameter, mm, $\pm 5\%$	315
Shielded power cable	Weight, kg, $\pm 1\%$	nineteen
	Wire length, mm, $\pm 1\%$	0.223
	Wire diameter, mm, $\pm 5\%$	3000
Key	Weight, kg, ± 0.001 kg	5.5
	Size, mm, ± 1 mm	0.007
Curved handle for RES electrode	Weight, kg, $\pm 1\%$	40x23x2
	Size, mm, $\pm 1\%$	0.253
	Wire length, mm, $\pm 1\%$	235 (length) x70 (width)
	Wire diameter, mm, $\pm 5\%$	2000
Curved handle for CAP electrode	Weight, kg, $\pm 1\%$	5.5
	Size, mm, $\pm 1\%$	0.253
	Wire length, mm, $\pm 1\%$	235 (length) x 70 (width)
	Wire diameter, mm, $\pm 5\%$	2000
Straight handle for CAP electrode	Weight, kg, $\pm 1\%$	5.5
	Size, mm, $\pm 1\%$	0.226
	Wire length, mm, $\pm 1\%$	210 (length) x 20 (width)
	Wire diameter, mm, $\pm 5\%$	2000
Remote Control	Weight, kg, ± 0.003 kg	5.5
	Size, mm, ± 1 mm	0,036
Remote control battery	Weight, kg, $\pm 10\%$	89x44x12
	Diameter, mm, $\pm 1\%$	0.003
Cylindrical Return Electrode	Weight, kg, ± 0.005 kg	20
	Size, mm, $\pm 1\%$	0,078
	Wire length, mm, $\pm 1\%$	121 (length) x 26 (diameter)
	Wire diameter, mm, $\pm 5\%$	2000
	The mass of the cable electrode return cylindrical, $\pm 0,005$ kg	5.5
Electrode bag	Weight, kg, $\pm 1\%$	0,078
	Sizes, mm, $\pm 1\%$	0.186
Flat return electrode	Weight, kg, $\pm 1\%$	305x75x75
	Sizes, mm, $\pm 1\%$	0.3
	Wire length, mm, $\pm 1\%$	221.8 (width) x 296 (height)
	Wire diameter, mm, $\pm 5\%$	2000
RES electrode with a diameter of 35 mm	Diameter, mm, $\pm 1\%$	5.5
	Diameter of the shaft connected to the handle, mm, $\pm 10\%$	35
	Weight, kg, $\pm 1\%$	6
RES electrode with a diameter of 50 mm	Diameter, mm, $\pm 1\%$	0,045
	Diameter of the shaft connected to the	50

	handle, mm, $\pm 10\%$	
	Weight, kg, $\pm 1\%$	6

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
RES electrode with a diameter of 65 mm	Diameter, mm, $\pm 1\%$	65
	Diameter of the shaft connected to the handle, mm, $\pm 10\%$	6
	Weight, kg, $\pm 1\%$	0.114
RES electrode with a diameter of 90 mm	Diameter, mm, $\pm 1\%$	90
	Diameter of the shaft connected to the handle, mm, $\pm 10\%$	6
	Weight, kg, $\pm 1\%$	0.209
30 mm diameter CAP electrode	Diameter, mm, $\pm 1\%$	thirty
	Diameter of the shaft connected to the handle, mm, $\pm 10\%$	3,5
	Weight, kg, $\pm 1\%$	0,032
40 mm diameter CAP electrode	Diameter, mm, $\pm 1\%$	40
	Diameter of the shaft connected to the handle, mm, $\pm 10\%$	3,5
	Weight, kg, $\pm 1\%$	0,047
55 mm diameter CAP electrode	Diameter, mm, $\pm 1\%$	55
	Diameter of the shaft connected to the handle, mm, $\pm 10\%$	3,5
	Weight, kg, $\pm 1\%$	0.051
65 mm diameter CAP electrode	Diameter, mm, $\pm 1\%$	65
	Diameter of the shaft connected to the handle, mm, $\pm 10\%$	3,5
	Weight, kg, $\pm 1\%$	0,065
CAP electrode with a diameter of 80 mm	Diameter, mm, $\pm 1\%$	80
	Diameter of the shaft connected to the handle, mm, $\pm 10\%$	3,5
	Weight, kg, $\pm 1\%$	0,085
Longitudinal electrode RES	Size, mm, $\pm 1\%$	100 (length) x 14 mm (diameter)
	Diameter of the shaft connected to the handle, mm, $\pm 10\%$	6
	Weight, kg, $\pm 1\%$	0.119
Longitudinal CAP electrode	Size, mm, $\pm 1\%$	100 (length) x 14 mm (diameter)
	Diameter of the shaft connected to the handle, mm, $\pm 10\%$	3,5
	Weight, kg, $\pm 1\%$	0,046
Massage nozzle	Diameter, mm, $\pm 1\%$	65
	Weight, kg, ± 0.003 kg	0,025
USB stick	Size, mm, ± 1 mm	56x7x19
	Weight, kg, $\pm 10\%$	0.011
Conductive Cream Anticellulite Cream	Size, mm, $\pm 1\%$	163 (height) x 97 (diameter)
	Volume, ml, ± 5 ml	1000
Conductive Cream Face Cream	Size, mm, $\pm 1\%$	163 (height) x 97 (diameter)
	Volume, ml, ± 5 ml	1000
Conductive Cream Care Cream RES electrode with a diameter of 65 mm	Size, mm, $\pm 1\%$	163 (height) x 97 (diameter)
	Volume, ml, ± 5 ml	1000

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acch

SPECIFICATIONS

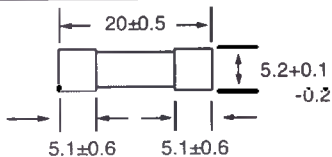No	Specifications		Value
ELITE NS, RGN			
1.	Input voltage, V		100-240 V
2.	Rated frequency		50/60 Hz
3.	Input fuse		2 * 4A, with a delayed response (T 4A L 250 V, in accordance with GOST IEC 60127-2-2013)
4.	Power consumption:		Standby time: 22 W ELITE NS - 360 VA (at 100 V), 352 VA (at 240 V) RGN - 150 VA (at 100 V), 145 VA (at 240 V) With an input voltage of 220-240V: ELITE NS -352 VA RGN - 145 VA
5.	Device protection from:		Random shorts between active and return electrodes Overheating
6.	Display		8 inches, 800x480 pixels with 10.1 " touch screen, 600 cd / m2
7.	Remote control battery		3 V, 220 mAh
8.	Connecting cable		3 pins per connector
9.	Shielded power cable		Input Connector Type - CEE7 / 7 Male Connector (Schuko) Output Connector Type - IEC 60320 C13 Connector
10.	USB stick		8 GB, USB 2.0
11.	The greatest effort required to engage the brakes on wheels, N		150
12.	The force required to move the cart rack with the main unit installed should not exceed, N		100
13.	Curved / straight handle for electrode CAP	The number of pins in the connector	5
		Connector type	Cylindrical
		Connector Diameter, mm, ± 1 mm	19
14.	Straight handle for RES electrode	The number of pins in the connector	4
		Connector type	Cylindrical
		Connector Diameter, mm, ± 1 mm	19
15.	Cylindrical Electrode Return	The number of pins in the connector	3
		Connector type	Cylindrical
		Connector Diameter, mm, ± 1 mm	19
16.	Fuse	Type of	Ceramic Cylinder Fusible 
		Number of fuses	2
		Sizes of a safety lock, mm	
		Fuse rated voltage, V	250
		Rated current of a safety lock, And	4

		Fuse rated time, s	See picture below.
		Fuses can be replaced in case of blown fuses by the user: disconnect the device from the power supply, remove the blown fuse from the fuse compartment, and then insert a new one, close the fuse compartment cover.	

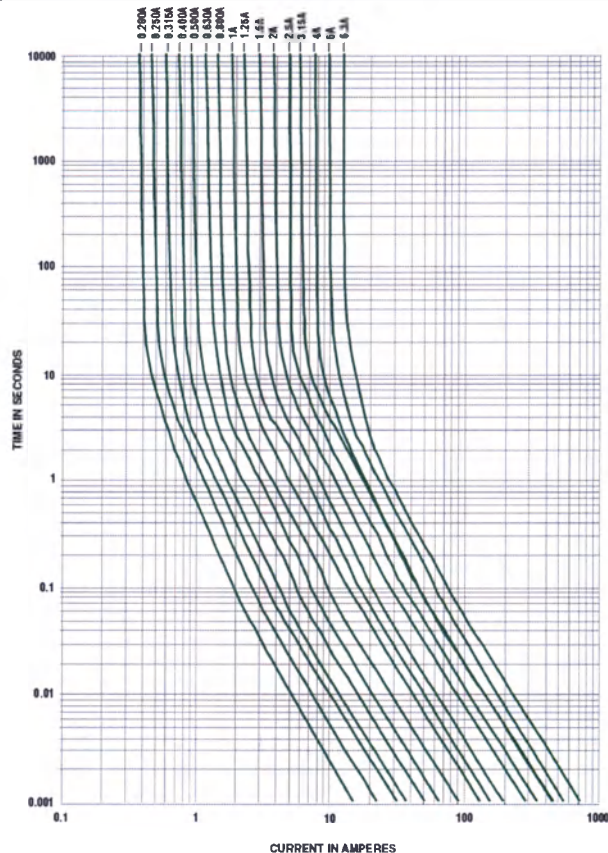
Fig. 69 Fuse response curve in seconds (vertical) versus current in amperes (horizontal).

ELITE NS		
	Resistive (RES) mode	Capacitive (CAP) mode
Output Waveform:	Sinusoidal, with and without modulation	
Type of modulation	Amplitude	
Total Harmonic Distortion (THD):	<2.5%	
Output frequency:	448 kHz	
Maximum output voltage:	120 V ~ (RMS)	620 V ~ (RMS)
Maximum power output:	200 VA	450 VA
Maximum Locked Output Power:	200 watts	150 watts
Permissible error of the output frequency:	± 1 kHz	
Permissible error of output power:	± 10% - ± 2 W, or the larger of the two values	

RGN		
	Резистивный (RES) режим	Емкостный (CAP) режим
Output Waveform:	Sinusoidal, with or without modulation	
Type of modulation	Amplitude	
Total Harmonic Distortion (THD):	<2.5%	
Output frequency:	448 kHz	
Maximum output voltage:	90 V ~ (RMS)	620 V ~ (RMS)
Maximum power output:	65 VA	250 VA
Maximum Locked Output Power:	65 watts	50 watts
Permissible error of the output frequency:	± 1 kHz	
Permissible error of output	± 10% - ± 2 Вт, либо большее из двух значений	

power:

The device is connected and disconnected from the power supply using the power switch on the rear panel. When the switch is in the "0" (Off) position, power is not supplied to the unit.

When checking or replacing fuses, disconnect the shielded power cable from the unit and turn the switch to the "0" position ("Off" position). Replaceable fuses must have the above specifications. This information can also be found on the markings located next to the fuse compartment.

Percentage of Conductive Creams

Conductive Cream Anticellulite Cream

Component	% content
Water	65,090
PEG-6 stearate	9,000
Vaseline oil	5,000
PEG-6 hydrogenated palm / palm kernel oil	5,000
Capric / Caprylic triglyceride	3,000
Glycol stearate	3,000
PEG-32 stearate	3,000
Propylene glycol	2,800
Sodium sulfate	1,200
Benzyl alcohol	0,800
Lanolin PEG-75	0,500
Dipropylene glycol	0,400
Perfumery fragrance	0,300
Myrtle arabic extract	0,200
Forskolia Coleus Root Extract	0,200
Horse Chestnut Extract	0,150
Triethanolamine	0,149
Salicylic acid	0,100
Glycerol	0,070
Sorbic acid	0,030
Fucus Bubble Extract	0,015

Conductive Cream Face Cream

Component	% content
Water	74,12180
PEG-6 stearate	7,68000
Vaseline oil	5,00000
Capric / Caprylic triglyceride	3,00000
PEG-6 hydrogenated palm / palm kernel oil	3,00000
Glycol stearate	2,16000
PEG-32 stearate	2,16000
Sodium sulfate	1,20000
Glycerol	0,04000
Salicylic acid	0,12000
Lanolin PEG-75	0,50000
Triethanolamine	0,17820
Sorbic acid	0,04000
Benzyl alcohol	0,80000

Conductive Cream Care Cream

Component	% content
Water	68,980
PEG-6 stearate	9,000
Vaseline oil	5,000
PEG-6 hydrogenated palm / palm kernel oil	5,000
Capric / Caprylic triglyceride	3,000
Glycol stearate	3,000
PEG-32 stearate	3,000
Sodium sulfate	1,200
Lanolin PEG-75	0,500
Benzyl alcohol	0,800
Perfumery fragrance	0,200
Triethanolamine	0,119
Salicylic acid	0,100
Glycerin	0,070
Sorbic acid	0,030

Characteristics of electrically conductive creams

	Conductive Cream Care Cream	Conductive Cream Face Cream	Conductive Cream Anticellulite Cream
Density, g / cm ³	1	1	1
pH	5,55	6,35	5,50
Viscosity, MPa * s	90000-150000	90000-150000	90000-150000
Electrical resistivity, Ohm * m	2,75 – 5,11	1,1 – 3,6	3,07 – 5,43
Maximum allowable contact time with skin, min	60	60	60

DISPOSAL

Dispose of the product in accordance with SanPiN 2.1.7.2790-10 for class B waste (epidemiologically hazardous waste). After hardware-based disinfection methods using physical methods and changing the appearance of the waste, eliminating the possibility of their reuse, class B waste can be accumulated, temporarily stored, transported, disposed of and disposed of together with class A waste (epidemiologically safe waste similar in composition to solid household waste waste). The packaging of Class B medical disinfected waste must be labeled to indicate that the waste has been decontaminated.

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.
.
.
.
.
.

Guarantees

INDIBA warrants that the product is free from defects in materials and workmanship.

The warranty on the device is valid for 2 (two) years from the date of purchase by a professional user. The warranty applies to any owner of the device during the warranty period in respect of defects in materials and workmanship.

Warranty on return and active electrodes, handles for electrodes depends on their frequency of use and care. The guarantee for return electrodes, handles for electrodes is two years, and for electrodes (capacitive and resistive) - one year.

The product has only been tested using conductive creams manufactured by INDIBA, therefore the warranty is valid only when using these conductive creams.

The manufacturer, assembler, installer or licensed importer is responsible for the safety, reliability and proper functioning of the product if:

1. All assembly, adjustment and repair operations are carried out by an individual or a company authorized by the manufacturer.
2. The electrical system in the treatment room complies with applicable requirements.
3. The device is used in accordance with the approved instruction manual.

The warranty does not cover the following cases:

1. Any malfunction or malfunction of the device if the integrity of the protective seal has been violated.
2. Any other malfunctions or malfunction of the device, not associated with defects in materials and workmanship, as well as improper use of the product.
3. Repairs made by an unauthorized service center or repaired using parts from unapproved INDIBA manufacturers.
4. Damage caused by the work of any party other than the manufacturer, dealer or authorized service agent

In order to receive technical support from a manufacturer, dealer or authorized service agent during the warranty period, the owner must:

1. Send a letter (by e-mail cs@indiba.com or by mail) to INDIBA, the dealer or an authorized service agent with a description of the problem during the warranty period.
2. If you need to send the device for inspection or repair, the attached letter indicates the serial number of the device (located on the label on the back of the device), the full name of the owner, full address, phone number and invoice number (if applicable) together with the name dealer or authorized service agent from whom the machine was purchased.

SERIAL NUMBER

The product identification mark is attached to the back of the machine. The serial number format is as follows:
YY S XXXX MM

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[Signature]

Troubleshooting

PROBLEM	ДЕЙСТВИЯ
"The product does not heat up"	• Make sure that all elements of the product, and especially the return electrode, are connected correctly. • If the indicator next to the start / stop button on the control panel is off, the start button has not been pressed. • If the indicator next to the start / stop button of the procedure blinks and the timer for counting the procedure time has been stopped, the device waits for contact between the electrode and the patient. • If the indicator next to the start / stop button of the procedure does not blink, and the timer for counting the time of the procedure counts the time, the device has detected contact with the patient. • If the problem persists, contact INDIBA Technical Support.
Capacitive electrode causes discomfort (acute pain)	• Make sure that there are no cracks or small holes on the surface of the electrode. Otherwise, replace the electrode. • During the procedure, the entire surface of the electrode should come into contact with the skin. If only a certain part of the electrode is in contact with the skin, all current passes through this small area, which can cause discomfort. - If the treatment area is uneven, use a smaller electrode. • Make sure that enough conductive cream has been applied to the skin. • If the problem persists, contact INDIBA Technical Support.
The ready light does not turn on when the power switch is pressed.	• • Ensure that the shielded power cable is not damaged and that it is correctly connected to the power supply. • • • Make sure that the power outlet is working correctly (for testing, connect another product to the connector). • • • Check that the two input fuses have not tripped. Otherwise, replace them with new fuses with the same characteristics (see section "Technical data"). • • • If the problem persists, contact INDIBA Technical Support.

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The control panel does not work (the indicator does not turn on when the start / stop button is pressed).

The device does not always respond to the buttons of the control panel, even if the indicator of the control panel turns on when you press the start / stop button of the procedure.

- The battery is low. Replace her.
- If the control panel still does not work after replacing the battery (the indicator does not turn on when you press the start / stop button of the procedure), make sure that the battery is in the correct position.
- If the problem persists, contact INDIBA Technical Support.
 - Move the remote control closer to the unit, as interference may interfere (the range of the control panel is 4 to 5 meters).
 - Check the battery level on the settings screen. If the battery is low, replace them.
 - If the problem persists, contact INDIBA Technical Support.

ERROR MESSAGES

When the machine detects an abnormal situation, it terminates the procedure and displays an error message. Each non-standard situation is encoded using an error number. The device displays the error number, a brief description and any additional information available.

If most errors occur, the only option for the user is to acknowledge the error. This is necessary to protect the user, patient, and the device itself. Then the device automatically turns off.

Some errors allow the user to continue using the machine after confirming the error.

If an error occurs, write down all the information displayed on the screen and contact INDIBA Technical Support.

Glossary of technical terms

ELECTROMAGNETIC COMPATIBILITY (EMC)

The ability of an electronic product not to interfere with other products and not to be sensitive to interference from other products.

ELECTRICAL SAFETY

A number of conditions that a product must meet in order to ensure the safety of people during normal use.

FREQUENCY

The number of repetitions of one cycle in one second. This concept refers to the frequency of alternating current. Unit: Hertz (Hz).

ELECTRICAL VOLTAGE

The electrical potential of the product or electrical system. The difference of electric potentials between two points of the circuit. Unit of measurement: Volt (V).

ELECTRICITY

A physical variable that expresses the amount of electricity flowing through a conductor over a given unit of time. Unit Type: Amp (A).

ELECTRICAL RESISTANCE

Resistance to current flow provided by the conductor. Unit: Ohm (Ohm).

ELECTRICAL IMPEDANCE

The relationship between the AC voltage applied to the circuit and the current generated by a load that does not have a pure resistance, but instead contains a capacitor. Unit: Ohm (Ohm).

CAPACITOR

A system of conductors separated by a dielectric (insulating) membrane that stores an electric charge. Under ideal conditions, it does not dissipate energy, but simply stores and then transfers it.

TOTAL POWER (VA)

The product of the voltage applied to the impedance and the current passing through it. Applicable when the load is impedance rather than ohmic resistance.

ACTIVE POWER (W)

The amount of energy consumed per unit of time that is dissipated into the patient's tissue.

ACTIVE ELECTRODE

The electrode in which the largest increase in temperature is observed due to current density.

ACTIVE METAL ELECTRODE (RESISTIVE)

Active metal electrode without an insulating layer.

ACTIVE METAL ELECTRODE (CAPACITIVE)

Active metal electrode with an insulating layer.

RETURN FLAT ELECTRODE

It is also an electrode, but it has a much larger surface area than the active electrode. Therefore, the current density in it is lower, so that the temperature in the surrounding tissue does not increase.

RETURN CYLINDER ELECTRODE

It is also an electrode, but it has a much larger surface area than the active electrode. Therefore, the current density in it is lower, so that the temperature in the surrounding tissue does not increase.

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EMC Information
Manufacturer's Manual and Declaration - Electromagnetic Emissions

The product is intended for use in an electromagnetic environment as defined below. The buyer or user of the product should ensure its use in the specified environment.		
EMC Tests	Conformity	Electromagnetic environment - instructions
Radio interference according to CISPR 11	Group 1	The product uses conductively coupled high-frequency energy necessary for the functioning of the product itself. Thus, the radiation of the product does not cause interference in neighboring electronic equipment.
Radio interference according to CISPR 11	Class B	The product is intended for use in premises for domestic purposes, as well as in rooms directly connected to low-voltage distribution electric networks (general purpose electric networks)
Harmonic components of current consumption in accordance with IEC 61000-3-2	Class A	
Compliance of voltage changes, voltage fluctuations and flicker caused by the vehicle with the standards of IEC 61000-3-3	Conformity	

The product is intended for use in an electromagnetic environment as defined below. The buyer or user of the product should ensure its use in the specified environment

The product is intended for use in an electromagnetic environment as defined below. The buyer or user of the product should ensure its use in the specified environment.			
Immunity Test	Test level according to IEC 60601	Compliance Level	Electromagnetic environment - instructions
Electrostatic discharges (ESD) according to IEC 61000-4-2	± 6 kV - contact discharge ± 8 kV - air discharge	± 6 kV - contact discharge ± 8 kV - air discharge	The floors of the room should be made of wood, concrete or ceramic tile. If the floors are covered with synthetic material, then the relative humidity should be at least 30%
Nanosecond impulse noise according to IEC 61000-4-4	± 2 kV - for power lines ± 1 kV - for input-output lines	± 2 kV - for power lines ± 1 kV - for input-output lines	The quality of electrical energy in the building's electrical network must be consistent with the typical conditions of a commercial or hospital setting.
High energy microsecond impulse noise according to IEC 61000-4-5	± 1 kV - when applying wire-to-wire interference ± 2 kV - when applying wire-to-ground interference	± 1 kV - when applying wire-to-wire interference ± 2 kV - when applying wire-to-ground interference	
Voltage dips, interruptions and voltage changes in input power lines according to IEC 61000-4-11	<5% UT (voltage dip) 95% UT) for 0.5 period; 40% UT (voltage dip 60% UT) for five periods; 70% UT (voltage dip 30% UT) for 25 periods; <5% UT (voltage dip) 95% UT) for 5 s	<5% UT (voltage dip) 95% UT) for 0.5 period; 40% UT (voltage dip 60% UT) for five periods; 70% UT (voltage dip 30% UT) for 25 periods; <5% UT (voltage dip) 95% UT) for 5 s	The quality of electrical energy in the building's electrical network must be consistent with typical commercial or hospital conditions. If the product user needs continuous operation under interruptions in the mains voltage, it is recommended that the product is powered by a battery or uninterruptible power supply.
Industrial frequency magnetic field according to IEC 61000-4-8	3 A / m	3 A / m	The magnetic field strength of the industrial frequency at the designated installation location should correspond to the typical conditions of a commercial or hospital setting.

Note: UT is the voltage level of the electric network until the test impact is applied.

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Immunity Test	Test level according to IEC 60601	Compliance Level	Electromagnetic environment - instructions
Conducted interference induced by radio frequency electromagnetic fields in accordance with IEC 61000-4-6	3 V (rms) in the band from 150 kHz to 80 MHz	3 V (rms value)	The distance between the portable / mobile radio communications used and any part of the product, including cables, must be at least the recommended spatial separation, which is calculated in accordance with the above expression as applied to the frequency of the transmitter. Recommended spacing: calculated by the formula: $d = 1,2\sqrt{P}$ for the frequency band from 80 kHz to 800 MHz: $d = 1,2\sqrt{P}$ for the frequency band from 800 kHz to 2.5 MHz: $d = 2,3\sqrt{P}$ where P is the nominal value of the maximum output power in W in accordance with the value set by the manufacturer; d - recommended spatial spacing, m The field strength during the propagation of radio waves from stationary radio transmitters according to the results of observations of the electromagnetic environment should be below the level of compliance in each frequency band. Interference may occur in the vicinity of equipment marked with the mark. 
Radiated radio frequency electromagnetic field according to IEC 61000-4-3	3 V / m in the band from 80 MHz to 2.5 GHz	3 V / m	

Note:

1. At frequencies of 80 and 800 MHz, a larger field strength value is used;
2. The above expressions are not applicable in all cases. The propagation of electromagnetic waves is affected by absorption or reflection from structures, objects and people;
3. Field strength during the propagation of radio waves from stationary radio transmitters, such as base stations of radio telephone networks (cellular / wireless) and terrestrial mobile radio stations, amateur radio stations, AM and FM broadcast transmitters, television transmitters, cannot be determined by calculation with sufficient accuracy. For this, practical measurements of the field strength should be carried out. If the measured values at the product location are higher than the applicable compliance levels, then the product should be monitored to verify its normal functioning. If during the observation process a deviation from normal functioning is detected, then additional measures must be taken, such as reorienting or moving the product.

Recommended separation values between portable and mobile RF communications equipment and the product

The product is intended for use in an electromagnetic environment in which the levels of radiated interference are controlled. The buyer or user of the product can avoid the influence of electromagnetic interference by providing a minimum spatial separation between portable and mobile radio frequency communications (transmitters) of the product, as recommended below, taking into account the maximum output power of communications			
Nominal maximum output power of the transmitter P, W	Spacing d, m, depending on the frequency of the transmitter		
	d = 1,2√P in the band from 150 kHz to 80 MHz	d = 1,2√P in the band from 80 to 800 MHz	d = 2,3√P in the band from 800 MHz to 2.5 GHz
0.01	0.12	0.12	0.23
0.1	0.37	0.37	0.73
1	1.20	1.20	2.30
10	3.70	3.70	7.30
100	12.00	12.00	23.00

Note:

1. At frequencies of 80 and 800 MHz, a larger field strength value is used;
2. The above expressions are not applicable in all cases. The propagation of electromagnetic waves is affected by absorption or reflection from structures, objects and people;
3. When determining the recommended values of the spatial separation d for transmitters with a nominal maximum output power not shown in the table, the nominal maximum output power P in watts specified in the transmitter manufacturer's documentation is substituted in the above expressions.

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Sección 2ª,

Content

05	1. System requirements
05	1.1 Functional warnings
06	2. INSTALLATION
07	2.1 First launch of software
09	2.2 Creating Users
11	3. OPERATION
11	3.1 Patients
12	3.2 Description of the patient (condition)
13	3.3 Procedure
16	3.4 History
18	4. MANUAL SETTING
19	5. PRINT
19	5.1 Insert custom logo
19	5.2 Printing a story
20	5.3 Status printing
20	5.4 Stages of printing
21	6. SETTINGS
21	6.1 User Management
21	6.2 Language
21	6.3 Version
22	7. DATA
22	7.1 Database Export
22	7.2 Importing a database
23	7.3 Database Integration
24	8. PROTOCOLS
24	8.1 Export Protocol
25	8.2 Import Protocol

IMPORTANT: Directive 95/46 / EC of the European Parliament and of the Council of 24 October 1995 guarantees the protection of individuals with regard to the processing of personal data and the free movement of such data. Compliance with the requirements of this directive is mandatory in the European Union. The country in which you are located, whether it is part of the European Union or not, may have approved its own legislation in this area. If so, these requirements must be observed.

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1. System requirements

A **monitoring program** has certain minimum system requirements.

Ensure that your environment meets these requirements before proceeding with the installation:

- Screen resolution of 1024 x 768 or higher for optimal display.
- Browser that supports HTML5 (Mozilla Firefox, Google Chrome, Safari).
- Computer with Wi-Fi connection. If your computer does not have a Wi-Fi module or the Wi-Fi module is not activated, contact the support service of the manufacturer of this computer.

Requirements for CPU, RAM, internal memory, etc. same as for the operating system. More internal memory may be required to store more patient data.

1.1 FUNCTIONAL REQUIREMENTS

- Ensure that the browser is not configured to periodically clear history data.
- Do not use private browsing in a browser. The database will not work.
- The browser must save history, data cache, and session variables.
- ◆ Always back up before deleting your browsing history.
- Each browser uses its own database.

1. Installation

The device also has a unique password for access to software for monitoring the implementation of procedures (hereinafter referred to as the monitoring program). This program should be downloaded to the specified address:

www.indibadeepcare.com/en/downloads/


and enter the required password, the monitoring program will be installed on the computer.

The program is downloaded in a compressed folder in .ZIP format. There are two files in the ZIP folder: one is the folder for copying, and the other is the same compressed folder that will be saved as a backup.

1. Drag the folder from the .ZIP file to the hard disk of the computer (it must be copied to the root directory of the memory).

2. In order not to waste time searching for a monitoring program, you can create a shortcut by going to the copy location of the folder, right-clicking and selecting "Create Shortcut".

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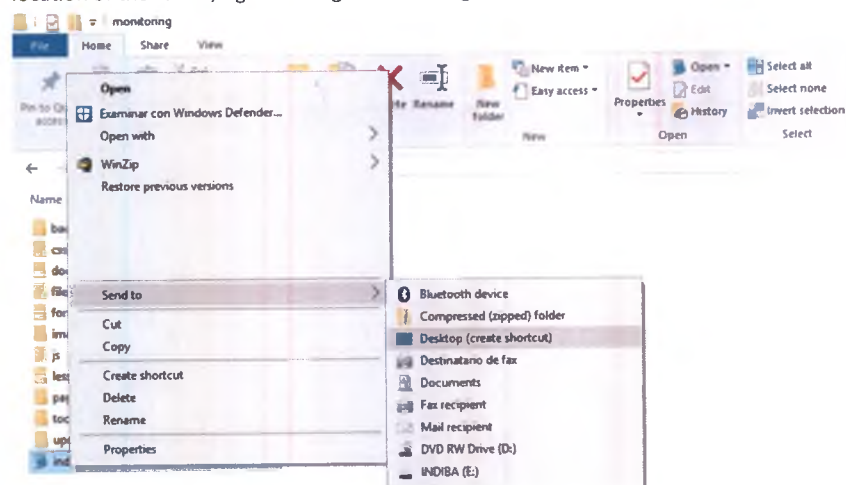
Handwritten signature


1.1 FIRST START-UP OF THE SOFTWARE

Before starting the monitoring program for the first time, it is necessary to establish a wireless connection between the computer and the device. To do this, follow these steps:

1. Click on the computer's Wi-Fi icon (it will flash white) and wait a few seconds.
 2. Connect the computer to the wireless network, which will be displayed as INDIBA_xxxxxxx, where xxxxxxx is the serial number of the machine.
 3. Enter a password to establish a connection. The password is wifiindiba.
 4. Once the connection is established, the Wi-Fi icon of the computer changes to blue and continues to flash.
- Now you can run the monitoring program on your computer for the first time. To do this, select the index.html file and run it. If several browsers are installed on your terminal (Chrome, Firefox or Safari), then "run as" (right-click on index.html) and "open with" appears, where you can select the desired browser. After a few seconds, the screen below will appear:


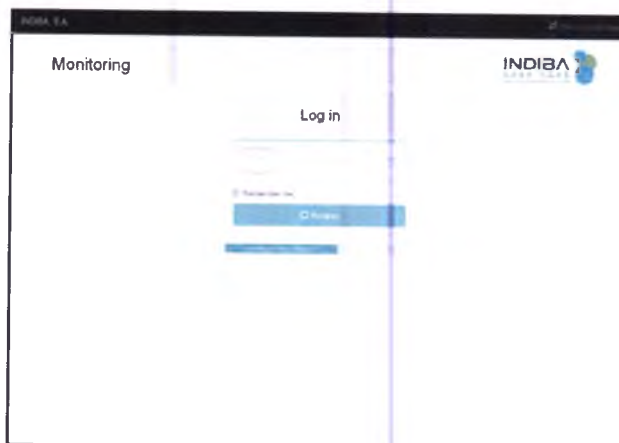
After confirming the creation of a database of ongoing procedures, the monitoring program will start and a window for entering the login and password by the user will appear.

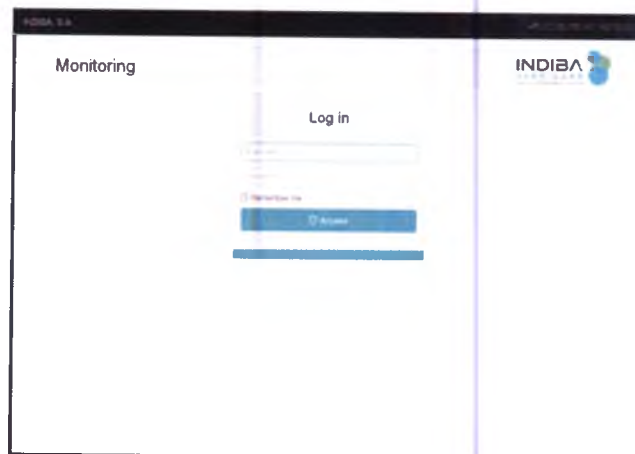At the first start, the program establishes a Wi-Fi connection with the device. If you have successfully connected wireless, the program expects to receive data from the device.

IMPORTANT: The user will not be able to work with the monitoring program until a connection is made to the device. After connecting, only a username and password are required to work with the program.

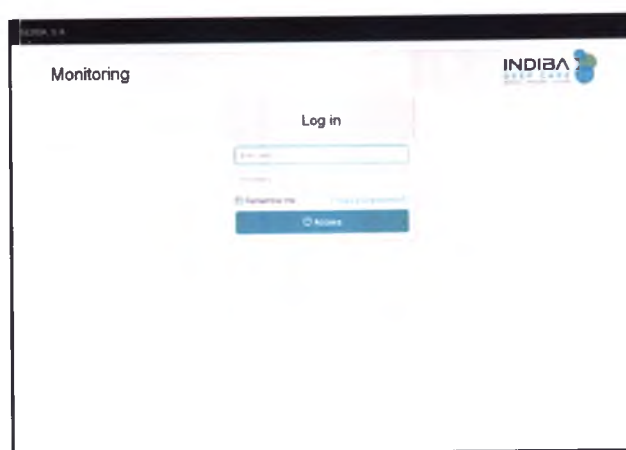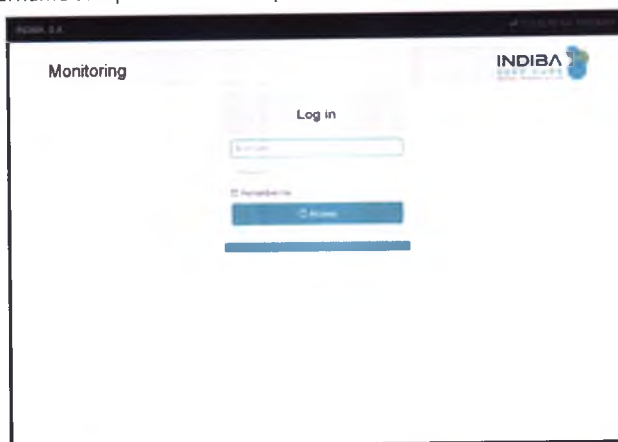The default username is indiba, the password is indiba.

IMPORTANT: if the “Remember me” checkbox is not selected, the user will need to enter their data each time to enter the program when it starts.

IMPORTANT: after entering with default data, the program detects that only this user exists, and instead of going to the patient window, redirects to the user creation window. This will happen every time a user enters the program until a new user profile is created instead of the default profile.

If the user has problems logging into the program with default data, the user can start working with the program by requesting a temporary password from the INDIBA customer service department. To do this, click on the link “Forgot your password?” And a pop-up window will appear containing the code. It must be reported to the INDIBA customer service department by calling the telephone number provided in the manufacturer’s contact information or by writing to cs@indiba.com, after which INDIBA specialists will send you a new generated password, valid for 7 days.

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

[Handwritten signature]

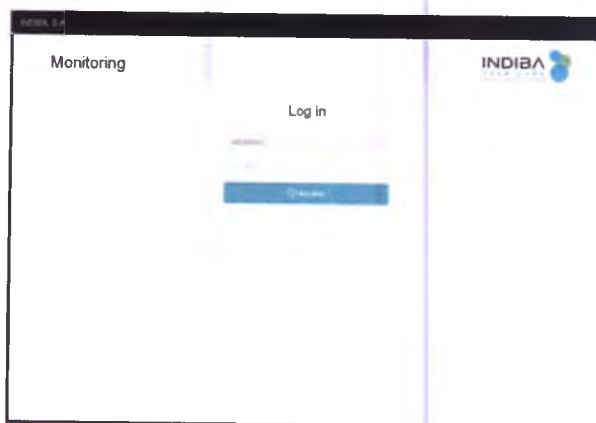
After entering the login and password, the program will open a user window.

1.1 CREATING USERS


The user creation screen appears immediately after entering the default username and password - indiba. This screen looks like this:

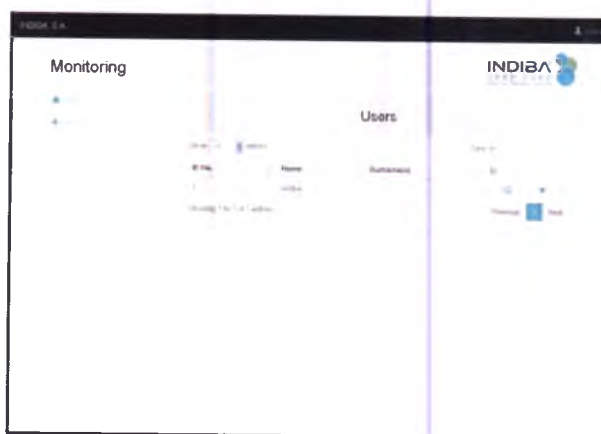The following explains how to create a new user named "John Simpson". To do this, fill in the fields of the pop-up menu, indicating the user identification number, first name, last name and password. Only lowercase letters are required to enter the password.

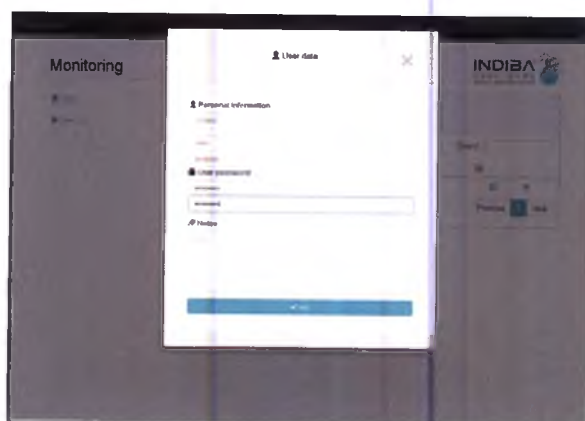After creating the user, a screen appears displaying all users in the database. On this screen you can edit a record or delete it. ✕

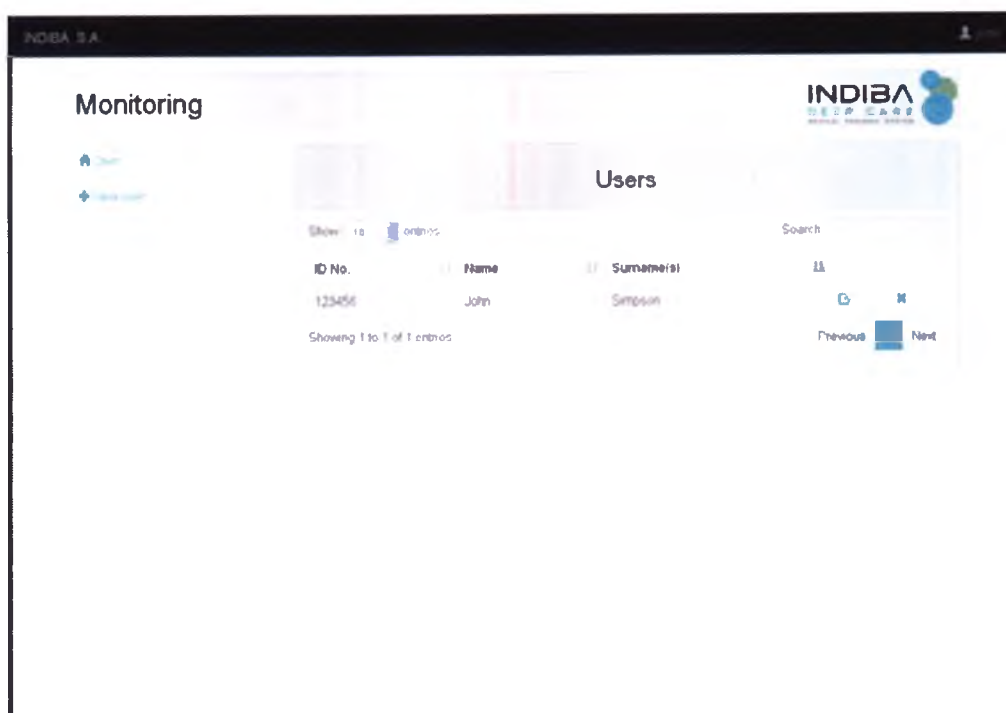If you do not want to create a new user, click on the “Start” link and go to the main screen (patients), using indiba data to enter the program. To create new users, see clause 7.1 “User management”, which describes this process in detail.

If you want to log in as a new user, you must log out of your Indiba account (in the black panel Logout”) and log in as a new user.


IMPORTANT: Each time you enter with standard indiba data, the program will redirect users to the screen, regardless of whether other users have already been created or not. If you log in as another user, the program will go to the main screen (section 3.1 Patients), where you can work with the program.

1. Operation

1.1 PATIENTS


The start screen is the patients screen. All operations with patients can be performed on this screen, in addition to procedures based on the amount of energy transferred and on the time of the procedure, for which only the data obtained on the device itself are used to create a report and results.

The start screen looks like this:

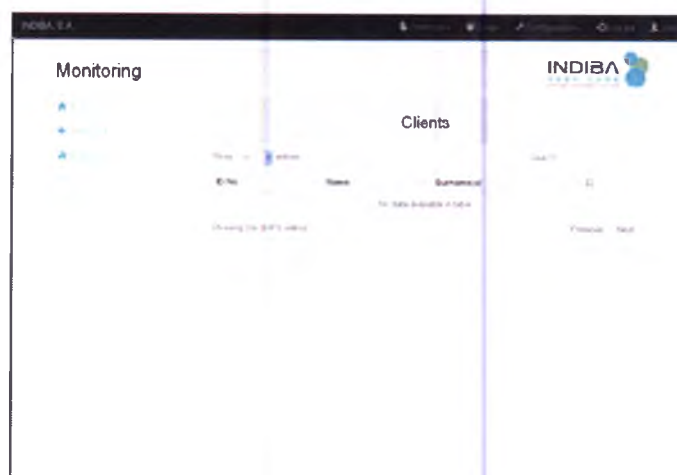At the top of the black bar is a drop-down menu ("Protocols" - PROTOCOLS, "Data" - DATA and "Configuration" - SETTINGS) and the option to log out of the current user's account ("Log out"). On the left, users can add new patients ("New Client") or perform the procedure without adding a new patient, but with data on the procedure ("Manual settings"). In the middle of the screen is a list of patients who can be searched directly and sorted in ascending or descending order by clicking the icons at the top.


To create a new patient, click "New Client" and the following pop-up menu will appear:

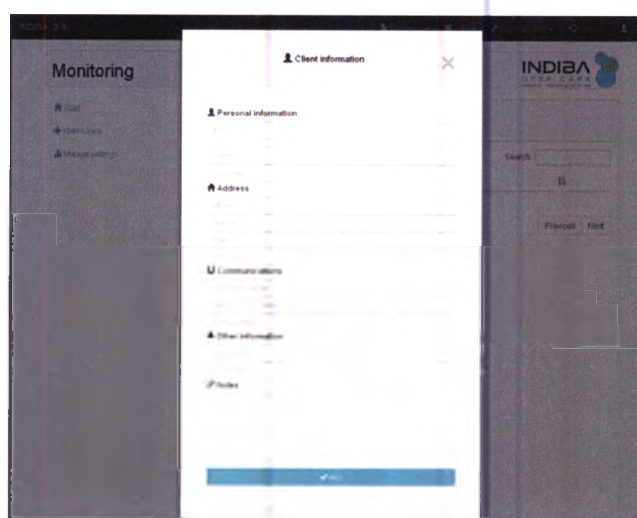Fill in the required fields and click "Add" - Add.

IMPORTANT: The [Identification Number] field cannot be changed after creation by the user. If the number is entered incorrectly, it must be deleted and re-entered. If you delete this user, all related data will be deleted.


After creating the patient, the following screen will appear:

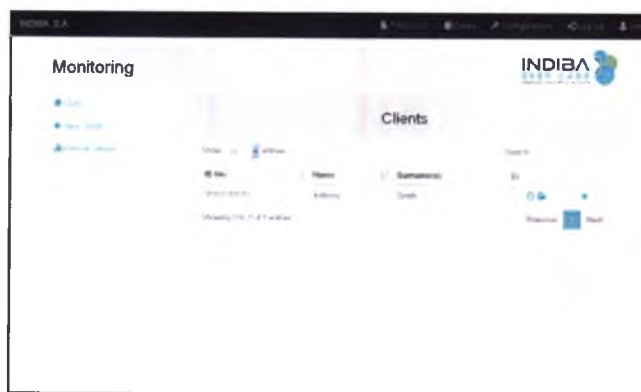
A series of icons will be displayed for each patient: the first to edit the record, the second to go to the patient's medical history and the third to delete the patient. After the patient is deleted, the data associated with it will be deleted from the database.

1.1 DESCRIPTION OF THE PATIENT (CONDITION)

Click on the name of the patient to view his description. The patient description screen has the same structure as the previous screen. However, the black toolbar at the top of the screen now only shows the username. Another change: the patient's name is displayed above the patient information table.

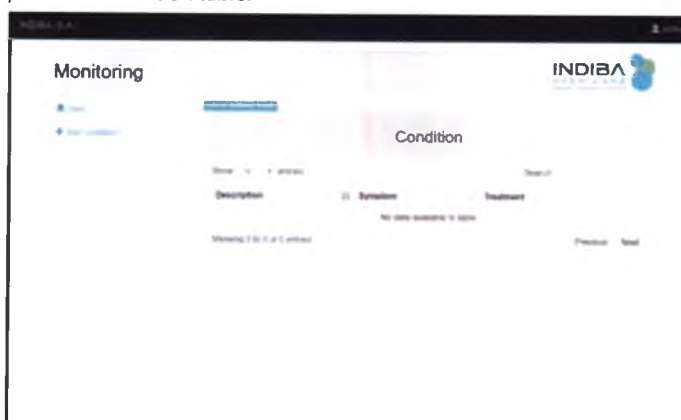
To add a description to the patient, click "Add condition" -Add a description, and a new pop-up menu will appear to enter the necessary information:

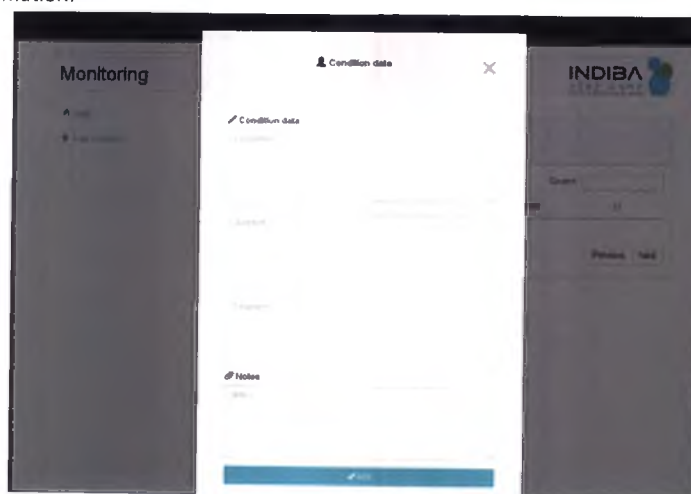Handwritten signature

After the information has been entered, click “Add” - Add”, and the following screen will be displayed:

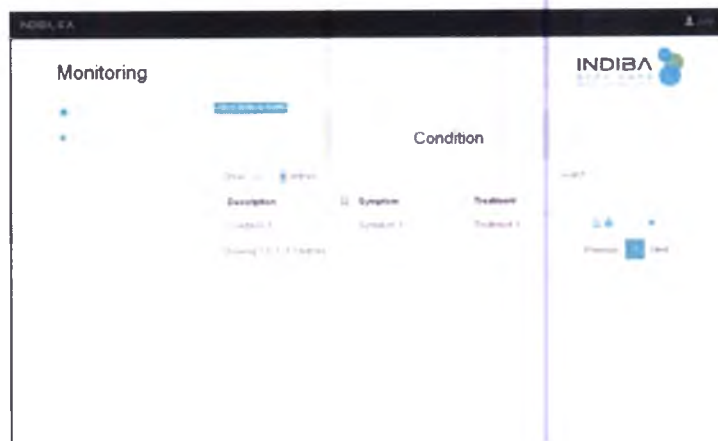
For each state, a series of icons will be displayed: the first to edit the record, the second to print the parameters of the procedure carried out for this state and the third to delete the state. When you delete a description, all related procedures are deleted from the database.


Click on a state to display the corresponding procedures that have been performed for that state.

1.1 PROCEDURE


The structure of the procedure screen is the same as the previous screens. The user can view the procedures that have already been completed or add new ones. The user name and corresponding status are displayed above the procedure table.

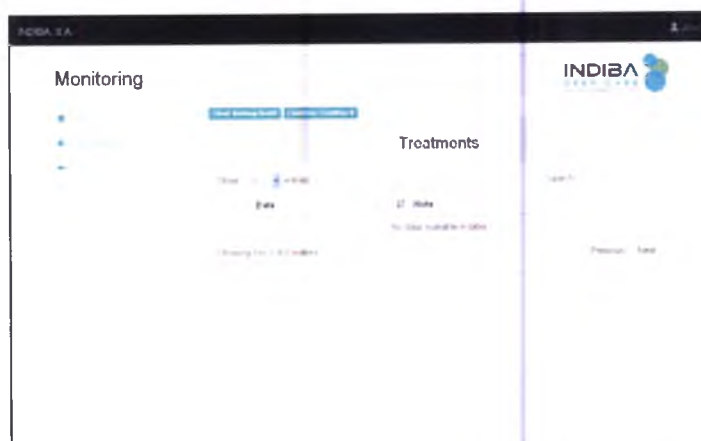
Click “Add treatment” - Add procedure to go to the next screen:

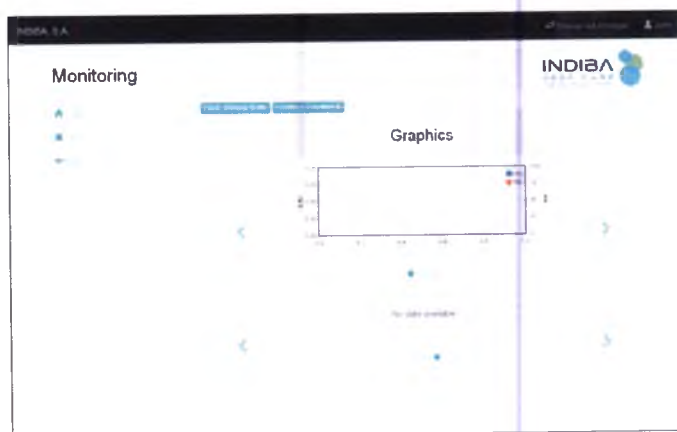At the top of the screen, in addition to the user name, a new message "Device out of range" - The device is out of range. If the device has a Wi-Fi connection with the computer, after a few seconds the name of the device and its serial number will appear at the top of the screen, as shown in the following figure:

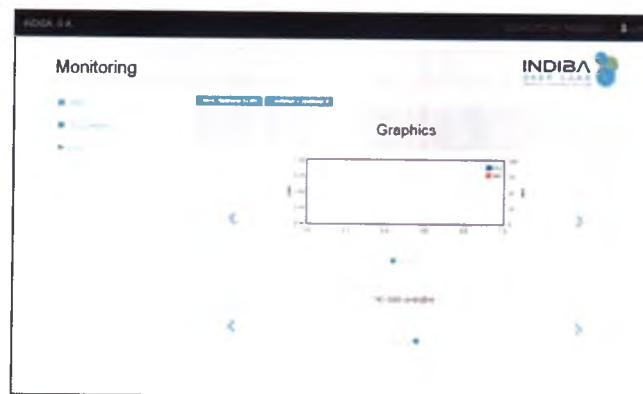
When the machine starts, the schedule changes over time, showing the last two minutes of the procedure. Click on the side arrows to display the data in a different format.

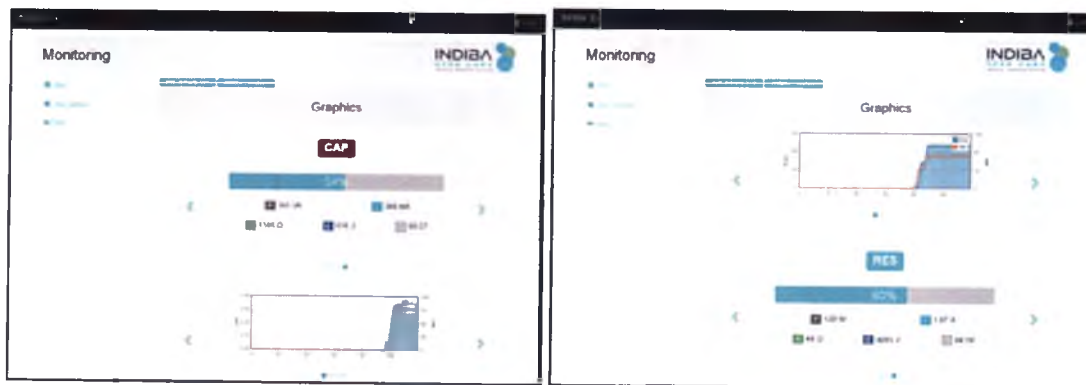
When the procedure moves to a new stage, a new schedule will be built for it. The switching of stages will continue until the user stops the procedure on the device or the data collection is stopped by pressing the "Stop capture" button - Stop data collection - in the program, in which case the data of the entire procedure will be displayed on the graph, and not just the last two minutes.

By clicking "Back" - Back, you will return to the previous page.


After receiving data about the procedure, the procedure screen will look like this:

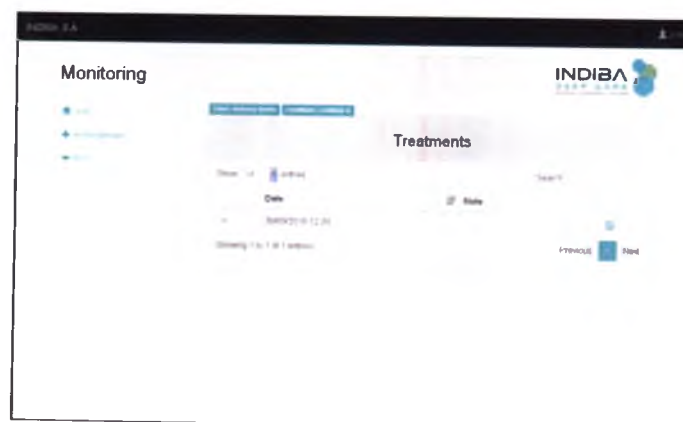
The Notes field remains blank at the end of the procedure, but you can edit it by clicking on the edit icon.

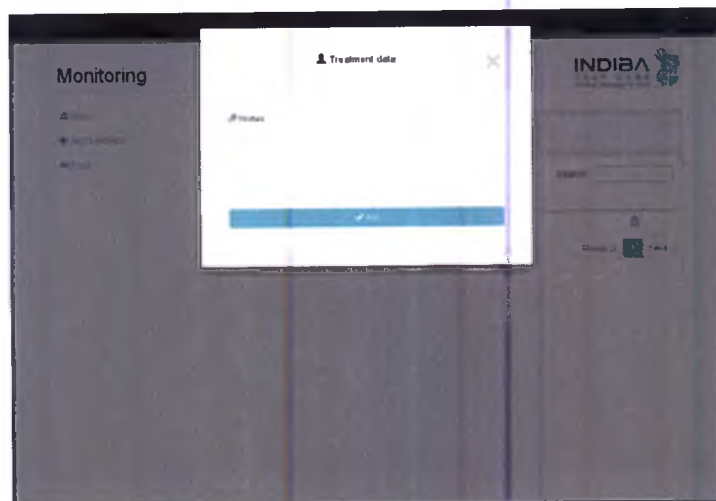On the procedures screen, there is a downward arrow that indicates that additional information can be revealed. After clicking on the arrow, the following will be displayed:

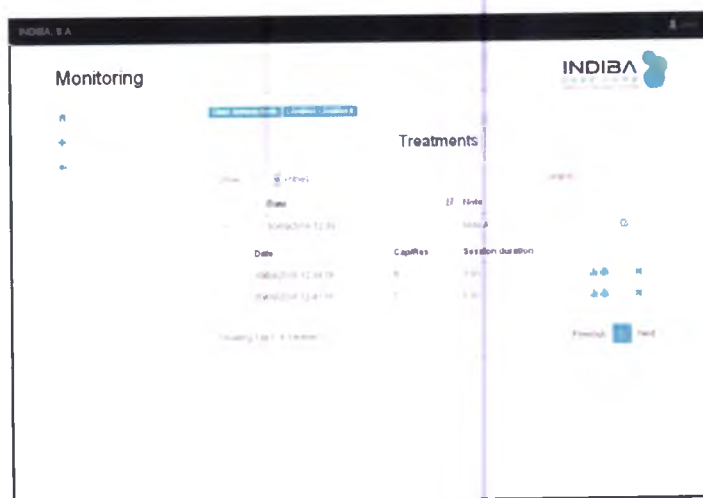
If the user clicks on a procedure and for some reason the data does not appear, or he clicks on a procedure that does not contain data, the following will be displayed on the procedure screen:

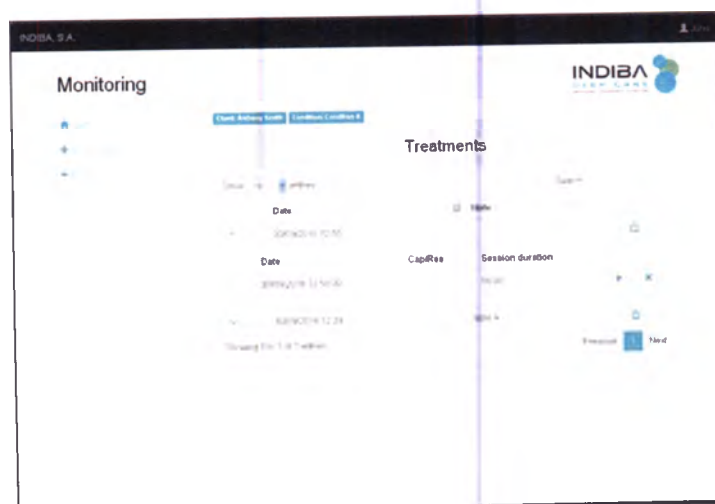
The procedure was interrupted. The duration of the session is 00:00. Click the icon to restart the session or to delete the session.

During the normal procedure, click the icon to display the results of this step, which can then be used for analysis.

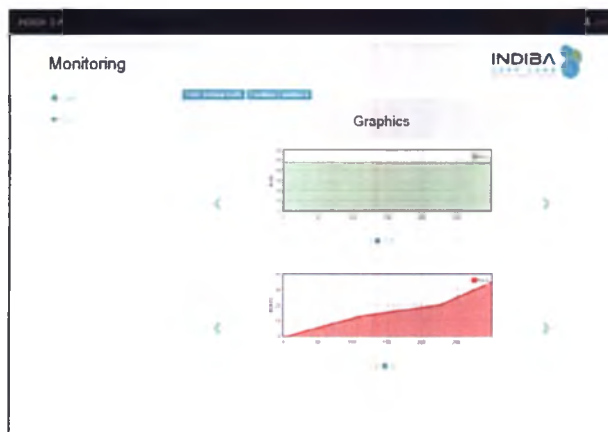
1.1 HISTORY

On the start screen, click on the icon next to the patient to access the history screen.

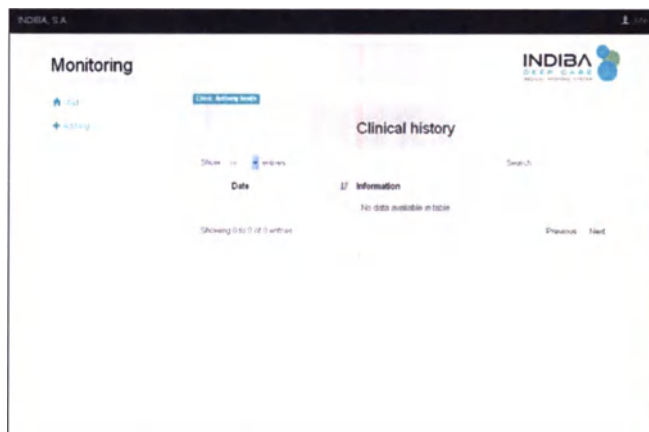
To add an entry, click "Add log" - Add a log, and a drop-down menu will appear to add the necessary information.

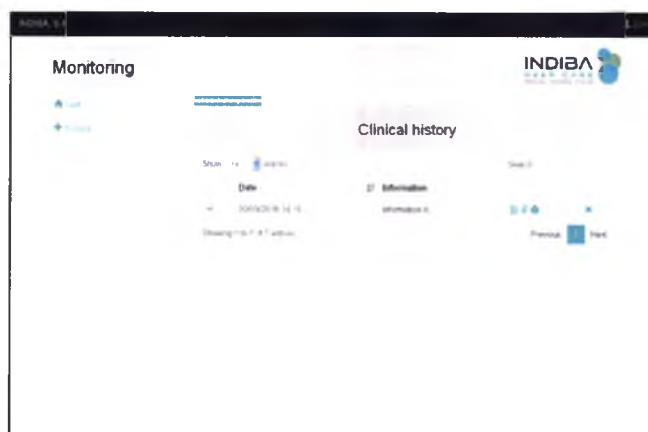
If the patient has data files, they can be linked to the record by clicking on the icon. All files provided by patients should be saved in a subfolder of the folder where the monitoring program was installed.


A screen appears that looks like this:

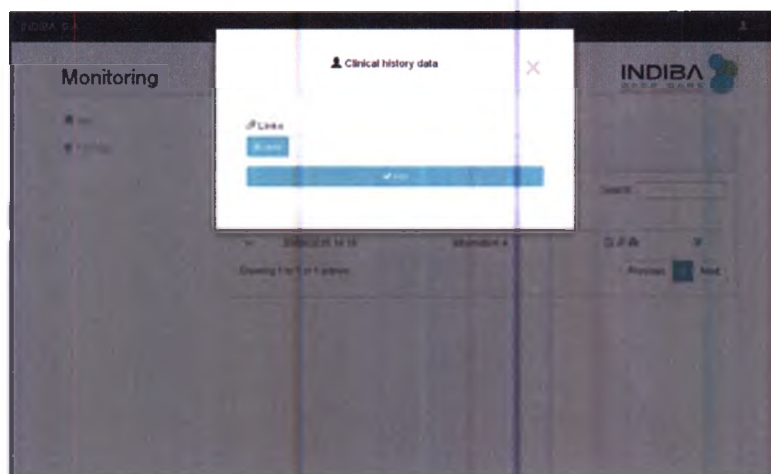Click Browse - Browse to select the desired file:

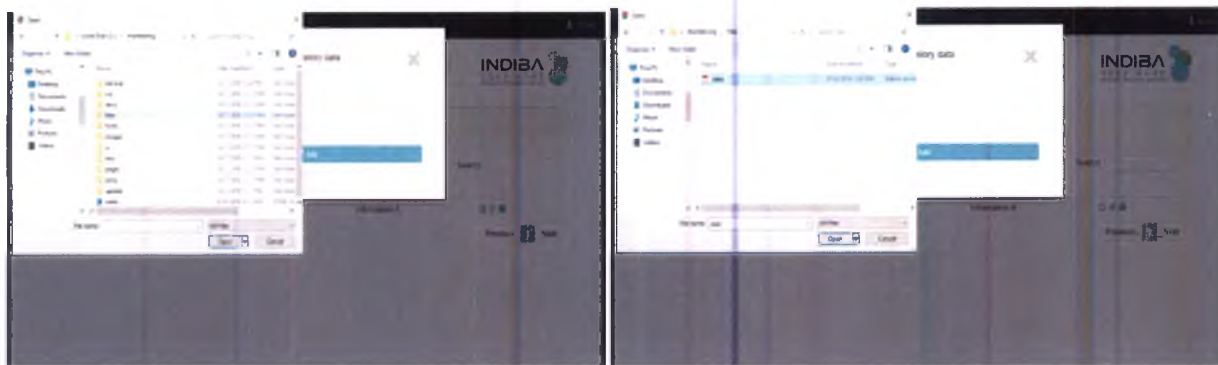After that, the entry will look like this:

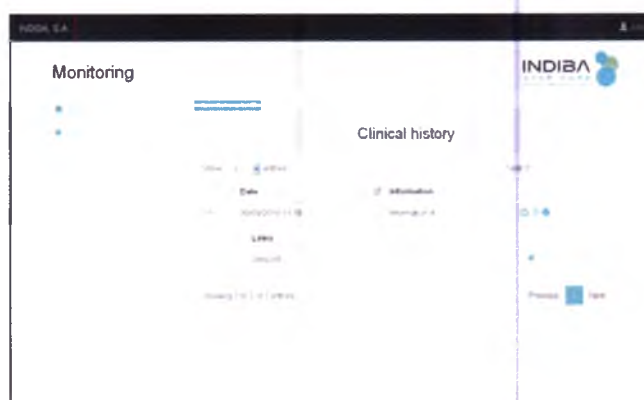
To delete an entry, all related files must be deleted first. If the required plug-in is installed in the browser, file links will open in a new window.

IMPORTANT: If a link to a file is deleted from the database, the linked file is NOT deleted from the computer. Files can only be permanently deleted in the folder in which they were saved.

1. Manual tuning

To perform the procedure for the patient without entering information about him into the program and provide a report on the procedure, click "Manual settings" - Manual setting. As shown in the following figures, the procedure is the same as the procedure for the saved patient.

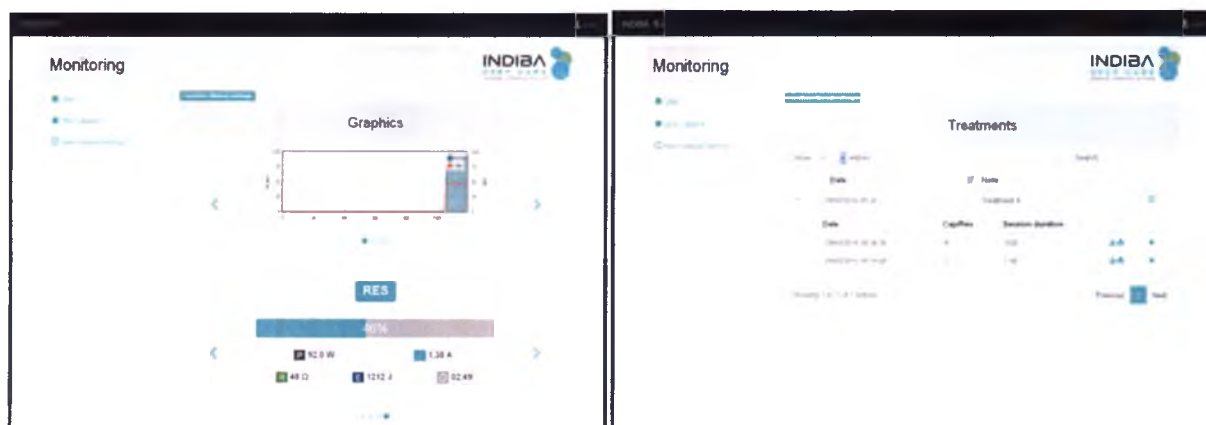
IMPORTANT: This procedure will not be saved in the shared database and will be lost when the user returns to the start screen or clicks "New manual settings" - New manual setting.

1. Printing

5.1 INSERTING A PERSONALIZED LOGO

The logo must be in jpg format and located in the “Images” folder (inside the “monitoring” folder) with the name logo.jpg. If the file height exceeds 100 points, the browser proportionally adapts the size.

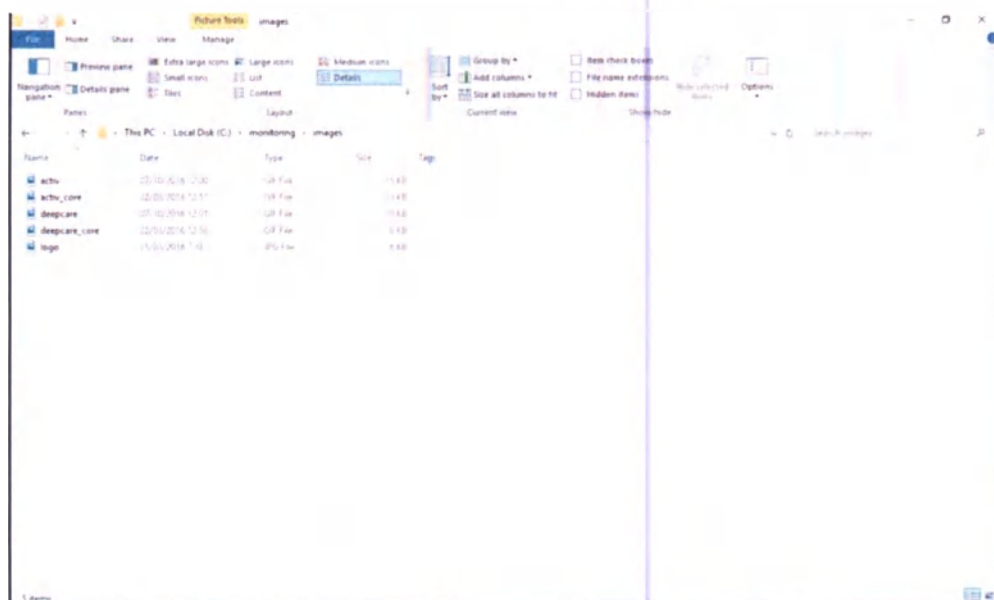The same image will be used in three reports prepared by the program.

5.2 PRINTING A HISTORY

Click on the icon in the desired history entry and the following screen will appear:


Click “Print” - Print, and a menu appears with a list of available printers.

IMPORTANT: The “Print” button will not be displayed on the printed document.


5.3 PRINTING STATE

To print the status, click "Print" on the status screen. The following screen will be displayed:

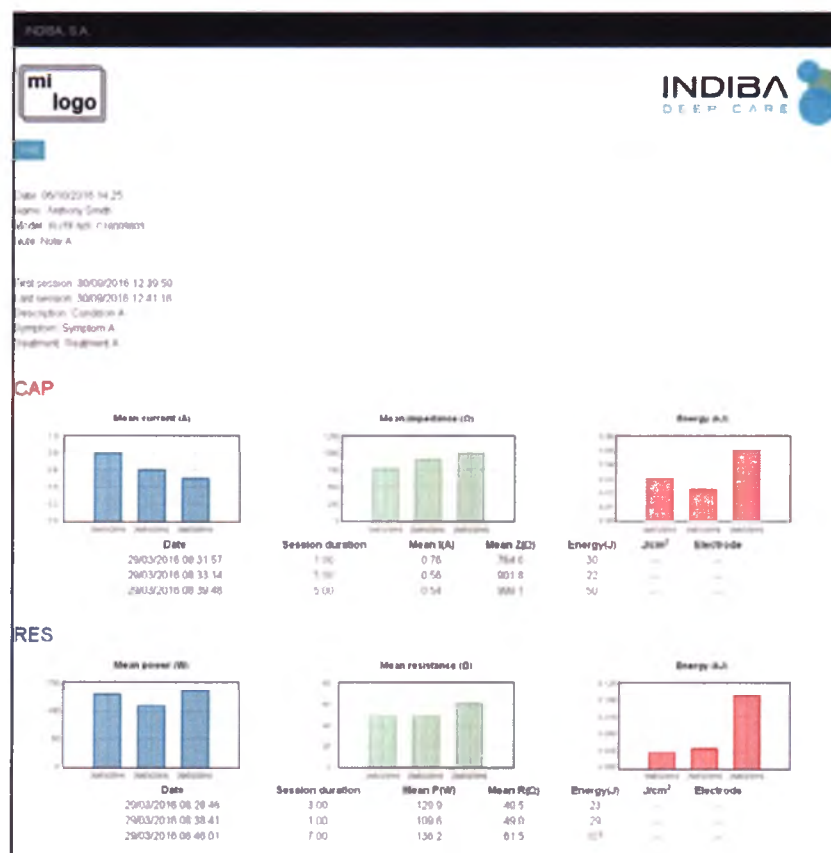5.3 PRINT STAGES


The following are sample reports for resistive and capacitive modes.

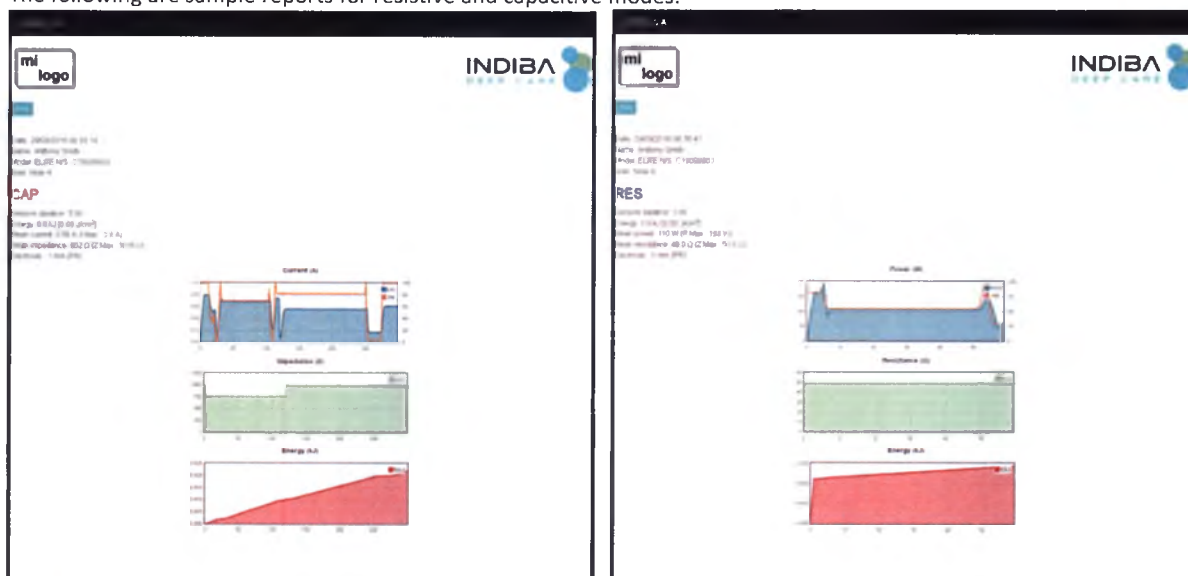7 Data

6.3 EXPORT DATABASE

Database backups should be done regularly to avoid any potential data loss. Warning message "ATTENTION! Remember to make a backup "will be displayed every fifteen days if a copy has not been made.

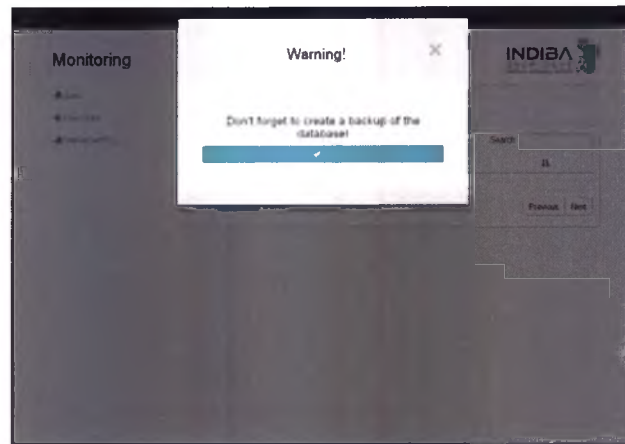IMPORTANT: The program displays only a message; automatic copying is not performed.


To make a copy, go to the data menu, click "Export DB" - Export DB (database), and a pop-up menu will appear.

Depending on the browser, the user will be asked to save the database in the program folder or use the default browser folder. It is recommended to save the file in the backup subfolder in the Monitoring folder so that it can be easily found if necessary.

IMPORTANT: The file name indicates the date the copy was created in the DDMMYY format.

6.3 DATABASE IMPORT

In order to export the database from another device, or restore the data, you can import the database. To do this, click on the "Import DB" - Import Database button, and the following screen will appear:

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
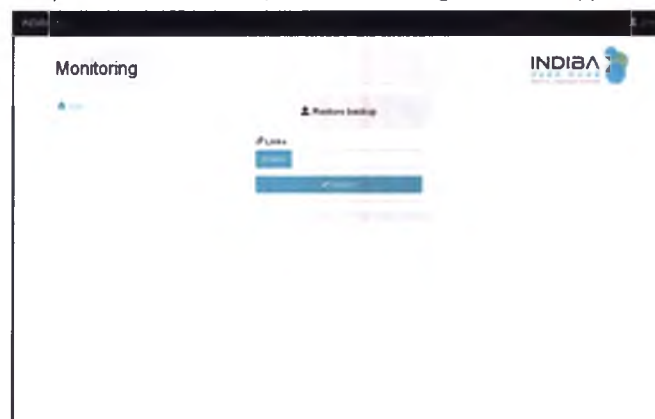
Handwritten signature

Click "Browse" - Browse to select a file to import the database. If the recommendations were followed, the file will be located in the "backup" subfolder in the folder where the monitoring program was installed.

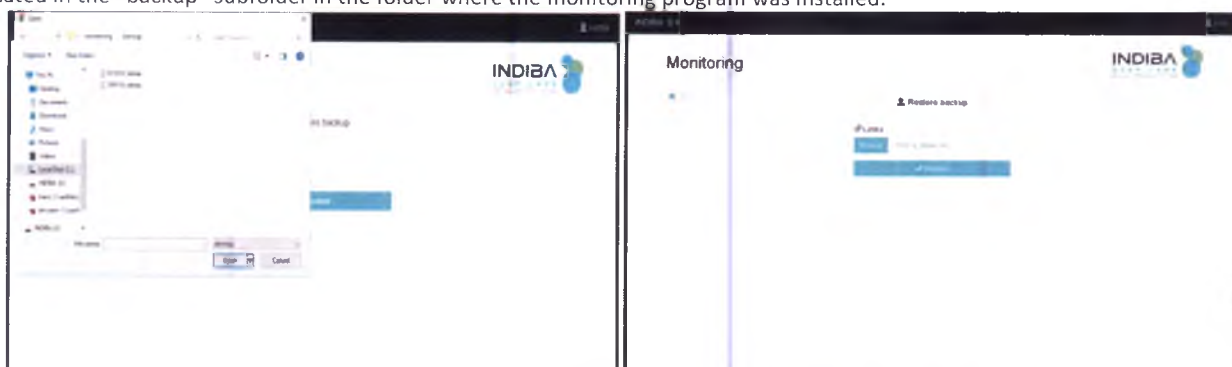
After the file is selected, you must click on the "Restore" button - Restore, the following confirmation screen will appear:

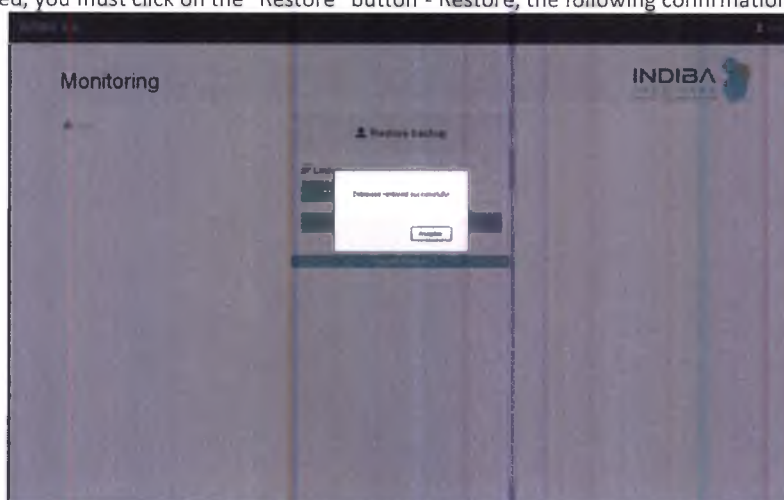IMPORTANT: after import, the new database will be overwritten on top of the old.

6.3 COMBINING DATABASE

This function was created for a situation when there are two monitoring programs on different devices. To combine these databases together on one device, the corresponding file must be exported, and then received and saved on another device. The subsequent steps are similar to those described in the previous section, and the user will receive confirmation on the screen if the process was successfully completed.

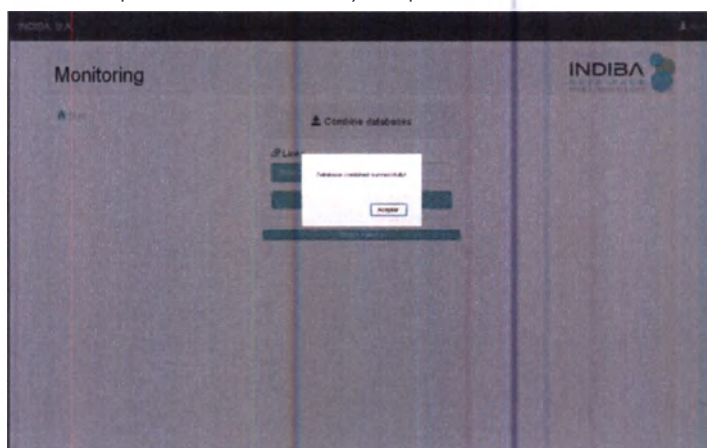
7 Protocols

The protocol menu allows you to import and export protocols created by the user on the computer, and vice versa.


6.3 EXPORT OF THE PROTOCOL

The tab "Export protocol" is used to send the user protocol from the computer to the machine. To do this, follow these steps:


1. Press the "Browse" button and select the folder for sending the program to the computer.

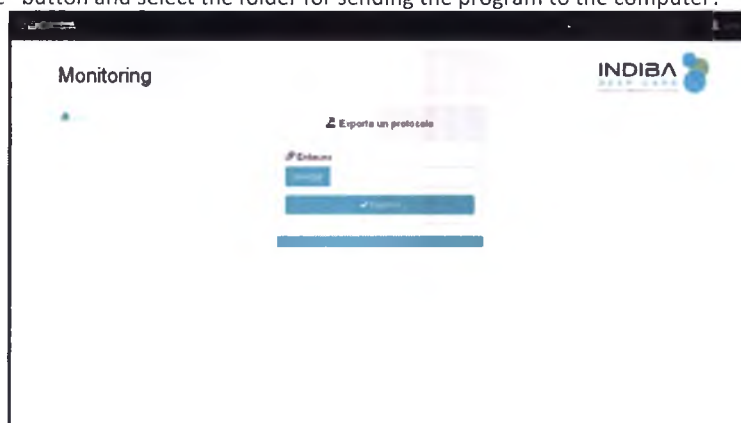2. On the device, open the page with user-saved procedures.

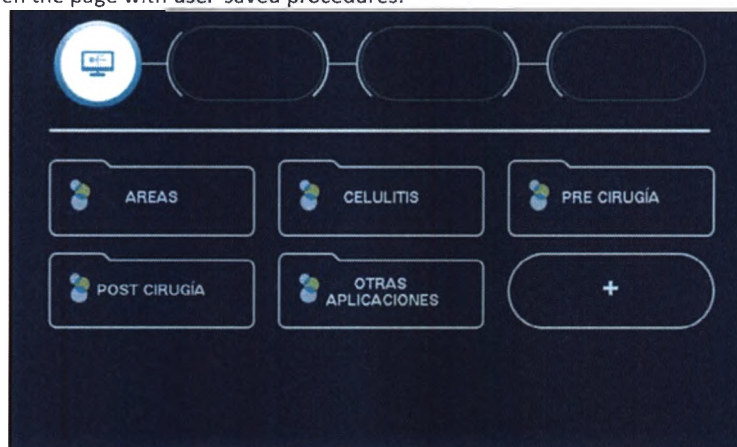
3. In the monitoring program, click the "Export" button. In a few seconds you will see a new program on the screen, its name is "TEST EXPORT".

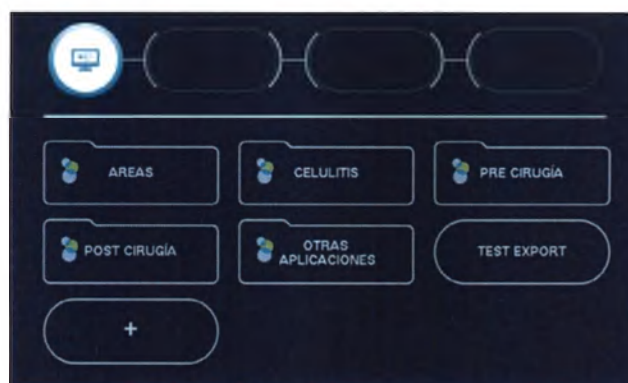
6.4 IMPORT OF THE PROTOCOL

The "Import protocol" tab is used to send the user protocol from the device to the computer. To do this, follow these steps:

1. On the device, select which program you want to transfer to the computer (in this case, a fake program called "TEST IMPORT" was used).

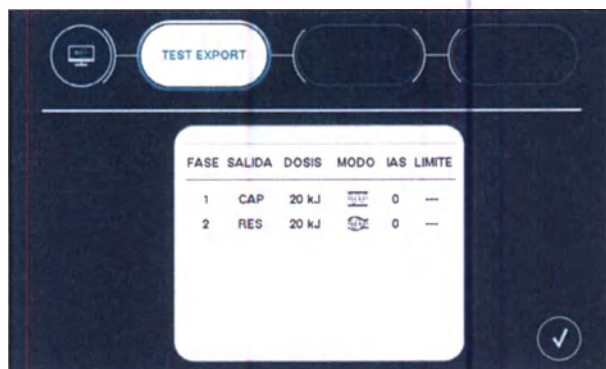2. In the monitoring program, click "Save" - Save.

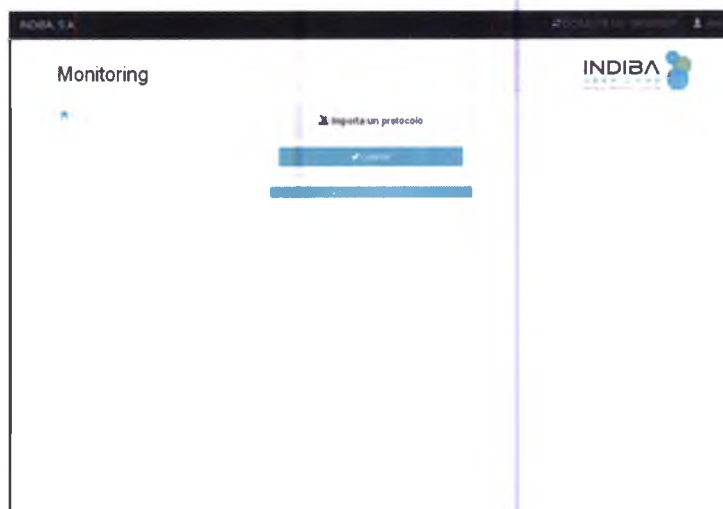3. Next, a pop-up window will appear to start the download. The file will be downloaded to the "Downloads" folder of your browser on your computer.

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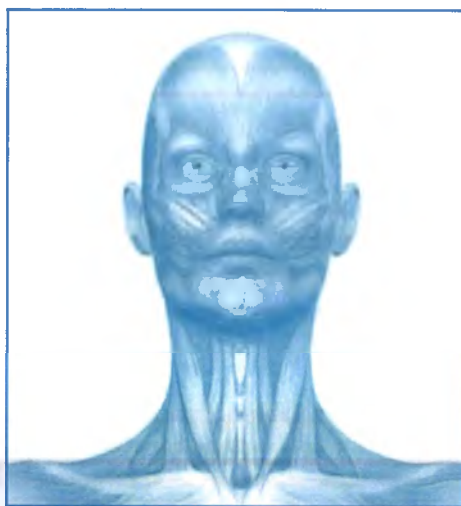
NECK AND FACE (30 min)

USED PRODUCTS

Conductive Cream Face Cream
Conductive Cream Care Cream
Face cream (not included)

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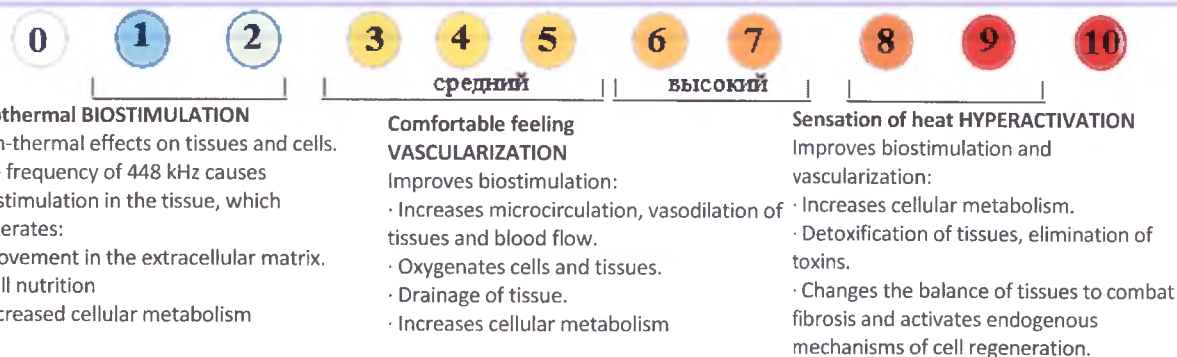
This protocol, which is exclusively owned by INDIBA, S.A., is protected by intellectual property law and may only be used by the rightful owners of the relevant equipment. Therefore, it is forbidden to reproduce, modify, distribute or transmit it in whole or in part, without the prior written consent of INDIBA, S.A. Все нарушения будут преследоваться по закону.

ARM (30 min)

Preparation of tissue for surgery, increased blood flow and oxygenation.

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Total: 4 sessions

NOTES

If other products are required, use them after the procedure.

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APPLICATION TECHNOLOGY

- 1. Position:** hold the electrode evenly on the surface of the body; do not tilt.
- 2. Pressure:** maintain sufficient pressure to ensure adequate contact between the skin and the surface of the electrode. The pressure should be stable and not cause discomfort to the patient.
- 3. Speed:** move the electrode in a circular motion, each time slowly increasing the diameter of the circle. Movements should be slow (one circle per second) and smooth.
- 4. Температура:** ensure that the temperature is suitable for the patient and is comfortable.
- 5. Гигиена:** To optimize the transfer of energy to the patient when using the product, it is recommended that the user hold onto the electrode handle during the procedure if it is necessary to use manual therapy methods, INDIBA representatives recommend the use of gloves (not included). This undoubtedly optimizes the transfer of energy to the patient and ensures the best result of the procedure.

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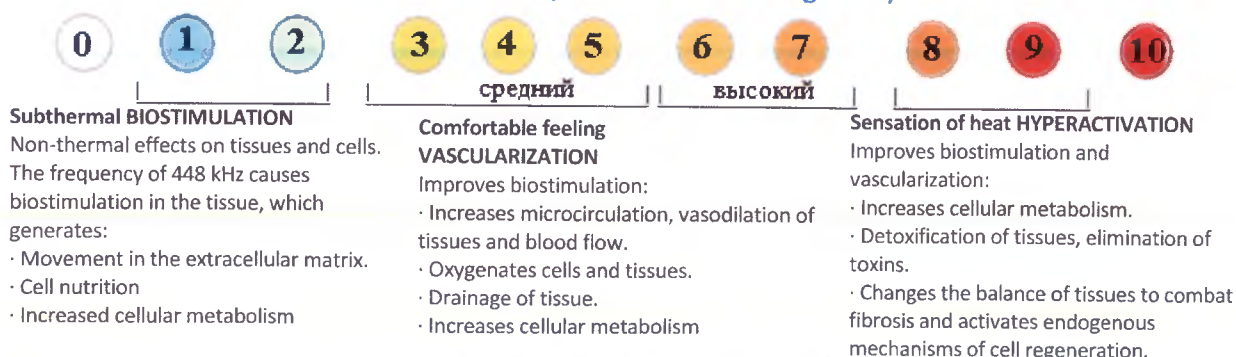
Preparation of tissue for surgery, increased blood flow and oxygenation.

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Total: 4 sessions

NOTES

If other products are required, use them after the procedure.

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97

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CHEST

INDIBA
DEEP CARE
MEDICAL FROGONIC SYSTEM

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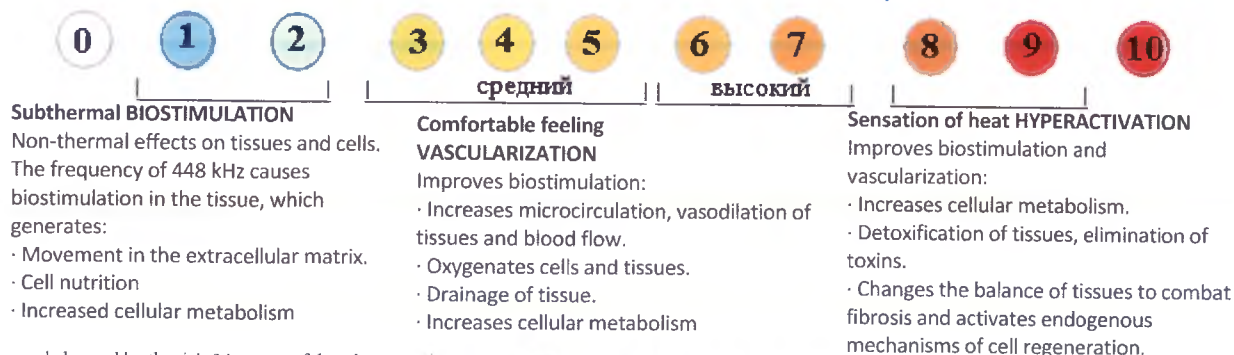
This protocol, which is exclusively owned by INDIBA, S.A., is protected by intellectual property law and may only be used by the rightful owners of the relevant equipment. Therefore, it is forbidden to reproduce, modify, distribute or transmit it in whole or in part without the prior written consent of INDIBA, S.A. измeнeнeй, рaспрoстpaнeнeй или пeрeдaвaтeль, eгo бeз пpeдвapитeльнoгo coглaсeния c INDIBA, S.A., выпyщeннoгo в письмeннoй фoрмe. Всe нapyшeния бyдyт прeдвoдeны пo зaкoнy.

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Preparation of tissue for surgery, increased blood flow and oxygenation.

POSITION

Lying on your back. Return flat electrode in the lumbosacral region.

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Total: 4 sessions

NOTES

If other products are required, use them after the procedure.

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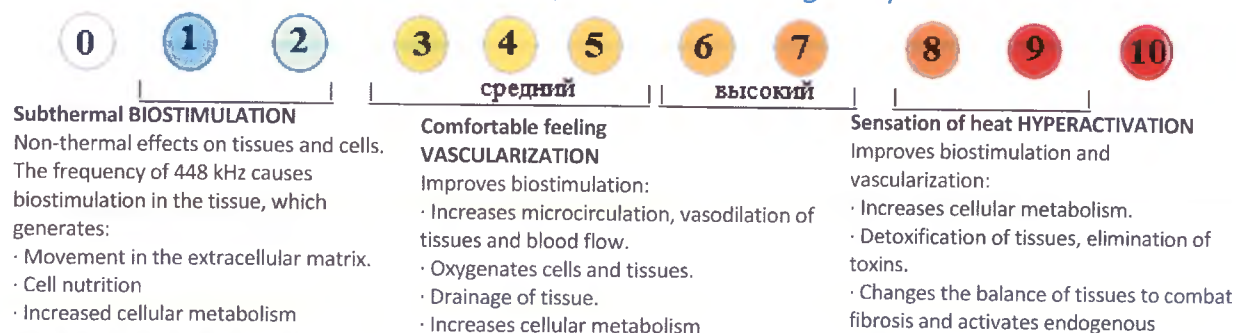
This protocol, which is exclusively owned by INDIBA, S.A., is protected by intellectual property law and may only be used by the rightful owners of the relevant equipment. Therefore, it is forbidden to reproduce, modify, distribute or transmit it in whole or in part without the prior written consent of INDIBA, S.A. Выпущено в письменной форме. Все нарушения будут преследоваться по закону.

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Preparation of tissue for surgery, increased blood flow and oxygenation.

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Total: 4 sessions

NOTES

If other products are required, use them after the procedure.

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Anti-inflammatory and decongestant treatment. Acceleration of the healing process. Tissue elasticity.

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CAP 2 min IAS* 1

APPLICATION DATA

Performing movements along the drainage line in the capacitive mode of the procedure. Ability to use at local drainage points.

SESSIONS

Once a week

Total: 4-5 sessions

NOTES

Treatment begins after suturing. In cases of sagging skin, it would be advisable to continue with the protocol from the group "Skin Tightening"

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
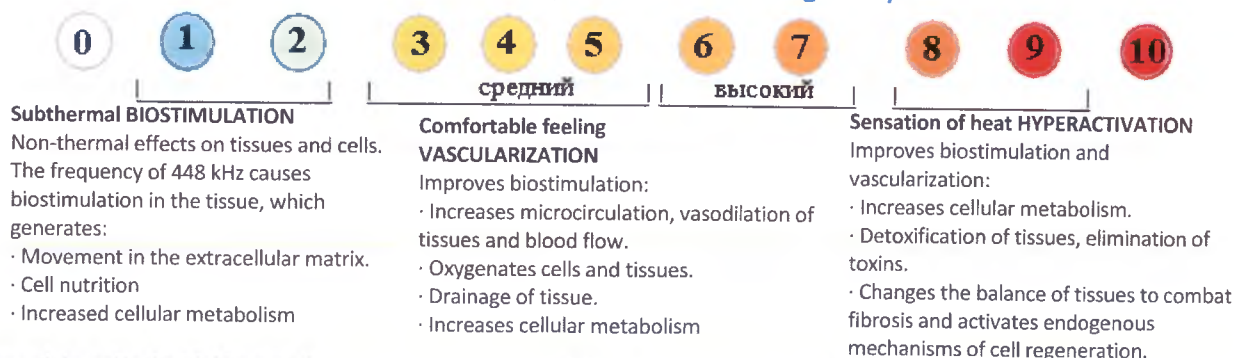
This protocol, which is exclusively owned by INDIBA, S.A., is protected by intellectual property law and may only be used by the rightful owners of the relevant equipment. Therefore, it is forbidden to reproduce, modify, distribute or transmit it in whole or in part, to change, disseminate or copy it without the prior written consent of INDIBA, S.A., issued in written form. All violations will be prosecuted by law.

BREAST IMPLANT INCAPSULATION_(30 min)

INDIBA
DEEP CARE
MEDICAL FRIGONIC SYSTEM

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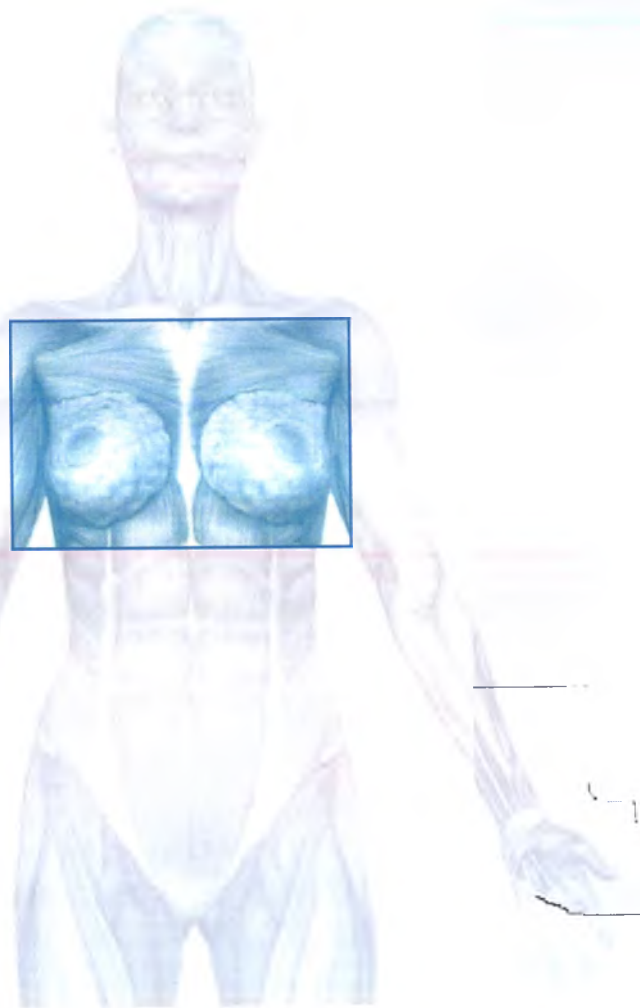SESSIONS

One or Two times a week

Total: 6 sessions to eliminate fibrosis

NOTES

Treatment begins after the appearance of periprosthetic fibrosis

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CAP 18 - 38 min **IAS*** 3 (depending on the size of the treatment area)

RES 2 min **IAS*** 3 (fixed on the sole of the foot in the treatment of feet)

APPLICATION DATA

Performing movements in the direction of drainage. It is possible to use at local drainage points.

SESSIONS

One-Two times a week

Total: 4-6 sessions

NOTES

During the procedure, it should be noted that hematomas can change color. Nodules after liposuction, palpable palpation, quickly normalize

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Treatment of sagging skin and fibrosis that remain after surgery.

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CAP 8 - 18 min IAS* 5 (depending on the size of the treatment area)

RES 10-20 min IAS* 7 (depending on the size of the treatment area)

RES 2 min IAS* 3 (fixed on the sole of the foot in the treatment of feet)

APPLICATION DATA

Performing movements in the direction of drainage. It is possible to use at local drainage points.

SESSIONS

Two times a week

Total: 10 sessions

NOTES

During the procedure, it should be noted that hematomas can change color. Nodules after liposuction, palpable palpation, quickly normalize

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Anti-inflammatory, decongestant and analgesic treatment.

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Lying on the abdomen: return flat electrode on the abdominal cavity

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CAP 10 - 15 min IAS* 3 (depending on the size of the scar)

CAP 5 min IAS* 1

APPLICATION DATA

Slow and smooth movements

SESSIONS

One or Two times a week

Total: 6 sessions

NOTES

If there are sutures, do not carry out the procedure directly on the scar, treat the area around

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APPLICATION TECHNOLOGY

1. **Position:** hold the electrode evenly on the surface of the body; do not tilt.
2. **Pressure:** maintain sufficient pressure to ensure adequate contact between the skin and the surface of the electrode. The pressure should be stable and not cause discomfort to the patient.
3. **Speed:** move the electrode in a circular motion, each time slowly increasing the diameter of the circle. Movements should be slow (one circle per second) and smooth.
4. **Температура:** ensure that the temperature is suitable for the patient and is comfortable.
5. **Гигиена:** To optimize the transfer of energy to the patient when using the product, it is recommended that the user hold onto the electrode handle during the procedure if it is necessary to use manual therapy methods, INDIBA representatives recommend the use of gloves (not included). This undoubtedly optimizes the transfer of energy to the patient and ensures the best result of the procedure.

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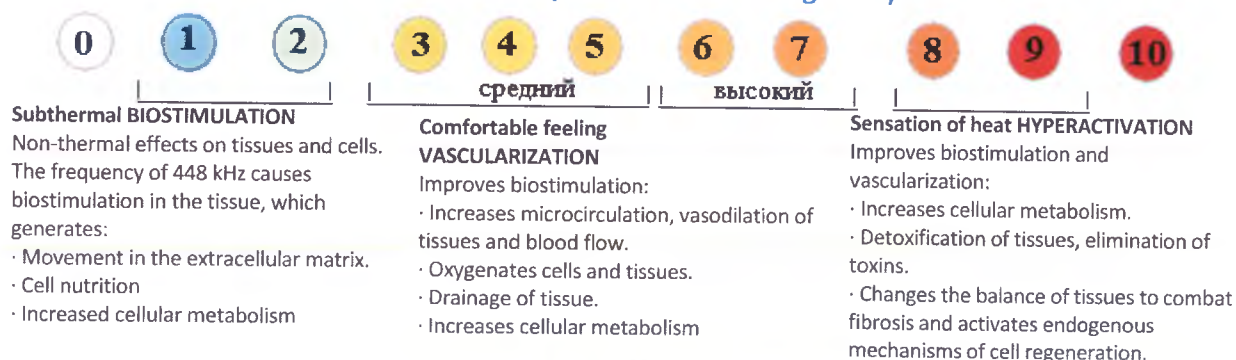
KELOID AND HYPERTROPHIC SCARS(10 - 20 min)

Anti-inflammatory, decongestant and analgesic treatment.

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Lying on the abdomen: return flat electrode on the abdominal cavity

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CAP 3 - 5 min IAS* 7 (depending on the size of the scar)

RES 5-13 min IAS* 7

CAP 2 min IAS* 3

APPLICATION DATA

Massage during a resistive treatment

SESSIONS

Once a week

Total:6 sessions

NOTES

Teach patient to do warm-up massage at home

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
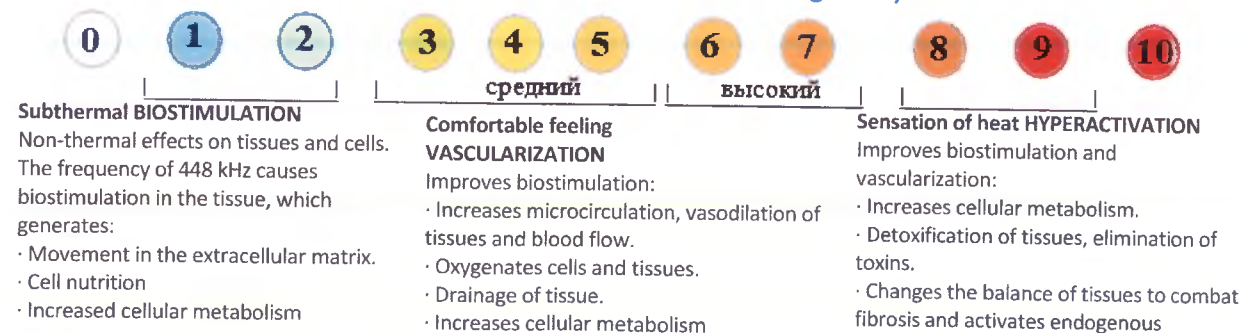
This protocol, which is exclusively owned by INDIBA, S.A., is protected by intellectual property law and may only be used by the rightful owners of the relevant equipment. Therefore, it is forbidden to reproduce, modify, distribute or transmit it in whole or in part without the prior written consent of INDIBA, S.A., issued in written form. All violations will be prosecuted by law.

BLEFAROPLASTY (20 min)

Decongestant, analgesic treatment to speed healing.

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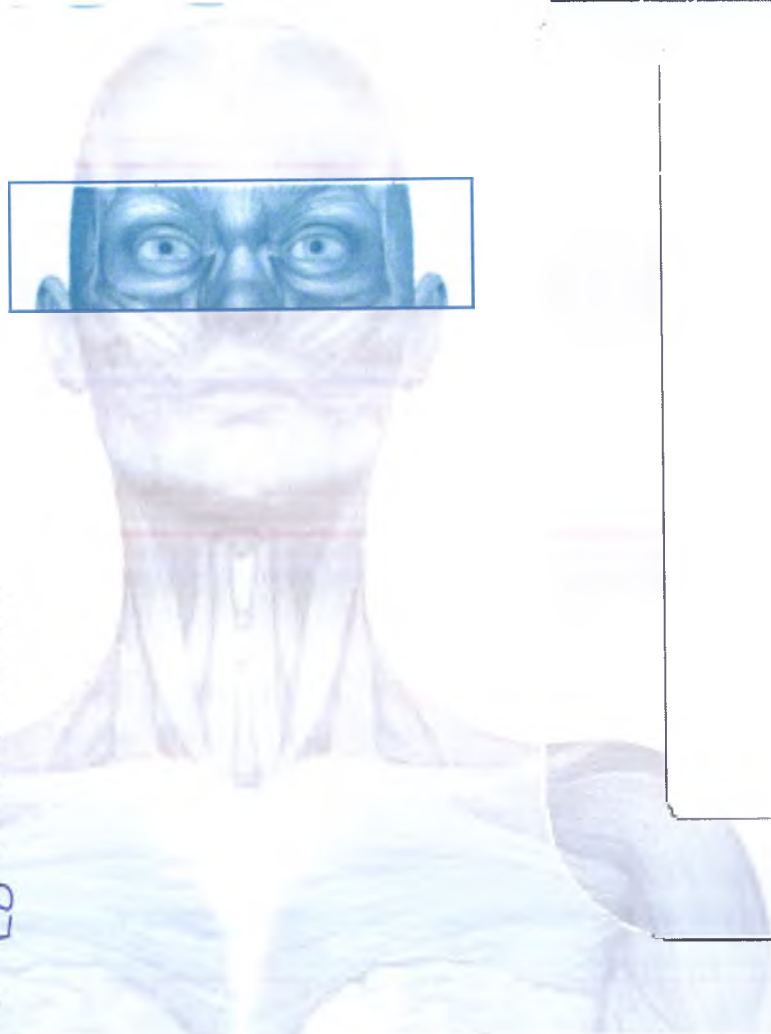CAP 8 min IAS* 3
RES 2 min IAS* 1 | EACH EYE

APPLICATION DATA

Slow, smooth movements in the direction of the temples and knots.

SESSIONS

Once a week
Total: 6 sessions

NOTES

If the RES electrode is uncomfortable for the patient, do not use it in initial sessions

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0	1	2	3	4	5	6	7	8	9	10
			средний			высокий				
Thermal BIOSTIMULATION non-thermal effects on tissues and cells. frequency of 448 kHz causes stimulation in the tissue, which operates: - Improvement in the extracellular matrix. - Full nutrition - Increased cellular metabolism			Comfortable feeling VASCULARIZATION Improves biostimulation: - Increases microcirculation, vasodilation of tissues and blood flow. - Oxygenates cells and tissues. - Drainage of tissue. - Increases cellular metabolism					Sensation of heat HYPERACTIVATION Improves biostimulation and vascularization: - Increases cellular metabolism. - Detoxification of tissues, elimination of toxins. - Changes the balance of tissues to combat fibrosis and activates endogenous mechanisms of cell regeneration.		

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Decongestant, analgesic treatment to speed healing.

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RES 5 min IAS* 1

APPLICATION DATA

Slow, smooth movements in the direction of the nodes.

SESSIONS

One-Two times a week

Total: 6 sessions

NOTES

If the RES electrode is uncomfortable for the patient, do not use it in initial sessions

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
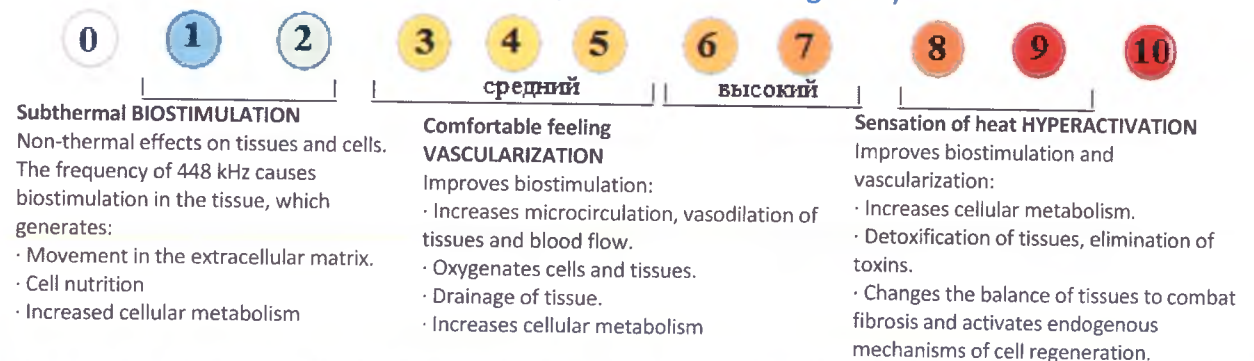
This protocol, which is exclusively owned by INDIBA, S.A., is protected by intellectual property law and may only be used by the rightful owners of the relevant equipment. Therefore, it is forbidden to reproduce, modify, distribute or transmit it in whole or in part without the written consent of INDIBA, S.A. Выпускено в письменной форме. Все нарушения будут преследоваться по закону.

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Decongestant, analgesic treatment to speed healing.

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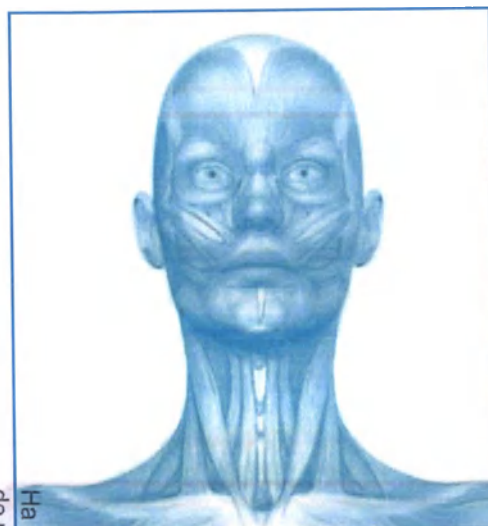Total: 6 sessions

NOTES

Treatment should begin three days after surgery, limiting the use of the chin, neck, face and eyelids. After carrying out the usual facial skin care procedures

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ABDOMINOPLASTY / LIPECTOMY (30 min)

Anesthetic and reabsorption treatment of hematomas and acceleration of the healing process.

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CAP 6 min IAS* 1

APPLICATION DATA

Drainage massage. Very slow and smooth movements

SESSIONS

Two times a week

Total: 10 sessions

NOTES

Check thermal sensations in the area that has not undergone surgery, conduct scar treatment in accordance with the existing protocol.

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APPLICATION TECHNOLOGY

1. **Position:** hold the electrode evenly on the surface of the body; do not tilt.
2. **Pressure:** maintain sufficient pressure to ensure adequate contact between the skin and the surface of the electrode. The pressure should be stable and not cause discomfort to the patient.
3. **Speed:** move the electrode in a circular motion, each time slowly increasing the diameter of the circle. Movements should be slow (one circle per second) and smooth.
4. **Температура:** ensure that the temperature is suitable for the patient and is comfortable.
5. **Гигиена:** To optimize the transfer of energy to the patient when using the product, it is recommended that the user hold onto the electrode handle during the procedure. If it is necessary to use manual therapy methods, INDIBA representatives recommend the use of gloves (not included). This undoubtedly optimizes the transfer of energy to the patient and ensures the best result of the procedure..

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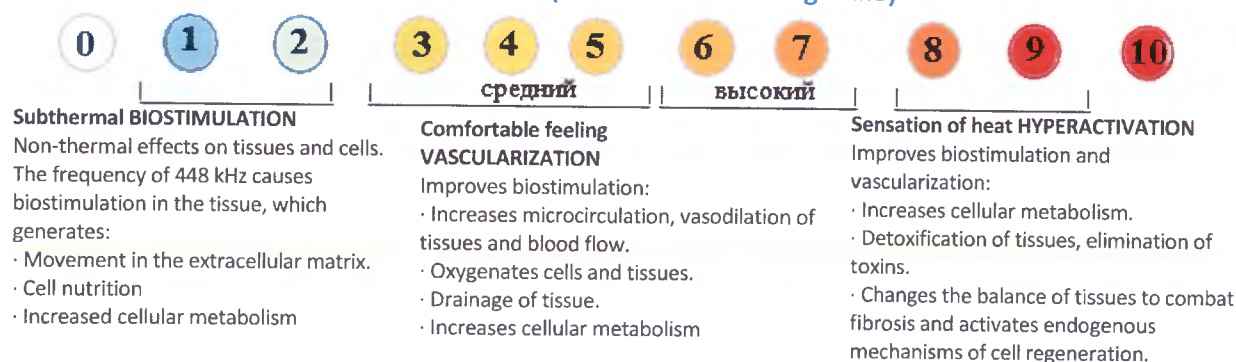
This protocol, which is exclusively owned by INDIBA, S.A., is protected by intellectual property law and may only be used by the rightful owners of the relevant equipment. Therefore, it is forbidden to reproduce, modify, distribute or transmit it in whole or in part without the written consent of INDIBA, S.A..

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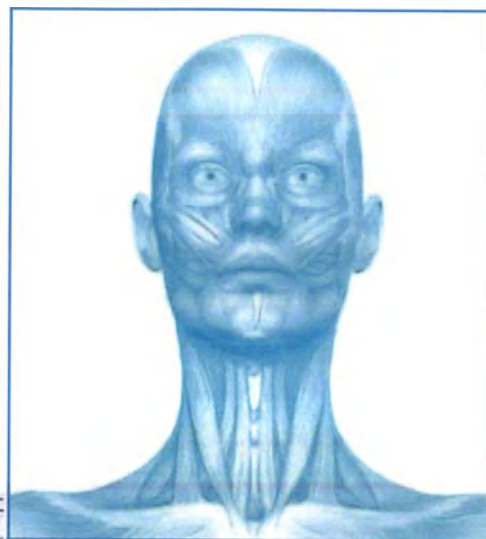
Conductive CreamCare Cream или
Conductive CreamFace Cream

RECOMMENDATIONS

It is recommended to carry out this procedure after targeted therapy or as concomitant therapy to improve the end result. Effectively combined with exercise and an appropriate diet.

COMMENTS

The proposed protocols are general and should be adapted to the individual requirements or circumstances of the patient, taking into account age, hormone imbalance, concomitant conditions, lifestyle, quality of life and history of surgical intervention. All of these factors can affect the end result. For external use only.


APPLICATION TECHNOLOGY

1. **Position:** hold the electrode evenly on the surface of the body; do not tilt.
2. **Pressure:** maintain sufficient pressure to ensure adequate contact between the skin and the surface of the electrode. The pressure should be stable and not cause discomfort to the patient.
3. **Speed:** move the electrode in a circular motion, each time slowly increasing the diameter of the circle. Movements should be slow (one circle per second) and smooth.
4. **Температура:** ensure that the temperature is suitable for the patient and is comfortable.
5. **Гигиена:** To optimize the transfer of energy to the patient when using the product, it is recommended that the user hold onto the electrode handle during the procedure. If it is necessary to use manual therapy methods, INDIBA representatives recommend the use of gloves (not included). This undoubtedly optimizes the transfer of energy to the patient and ensures the best result of the procedure.

CONTRAINDICATIONS

The procedure should not be carried out in the presence of a pacemaker or other electronic implants, in pregnant women, in case of skin wounds (open wounds or recent burns), thrombophlebitis and lack of sensitivity. Not contraindicated in the presence of implants

Patient sensations based on IAS value (IAS = INDIBA analog scale)

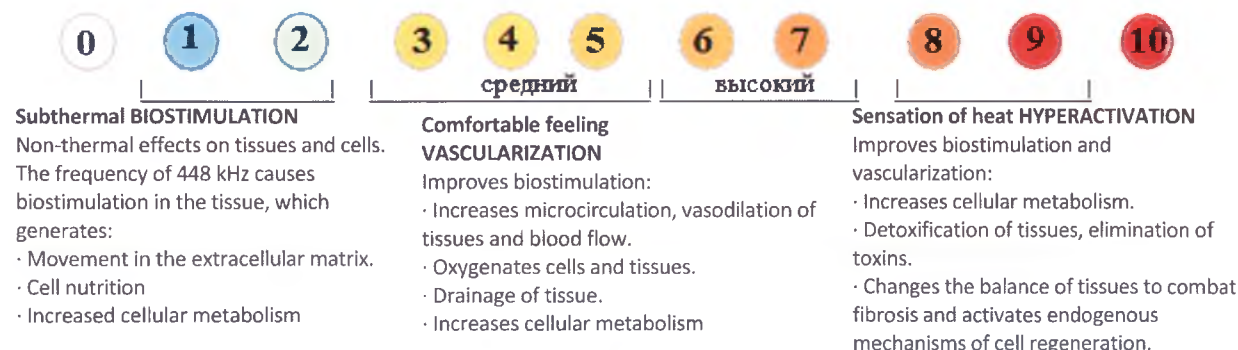NECK AND NECK LINE (30 min)

The procedure for increasing the elasticity and firmness of the skin. Epidermis regeneration and improved vascularization.

POSITION

Lying on your back. Return flat electrode between the blades.

DESCRIPTION OF STAGES OF THE PROCEDURE

CAP 12 min IAS* 7

RES 12 min IAS* 7

RES 6 min IAS* 5

APPLICATION DATA

Perform circular movements in the direction of the muscles and muscle cords. At the final stage (resistive mode) apply the methods of "tissue remodeling".

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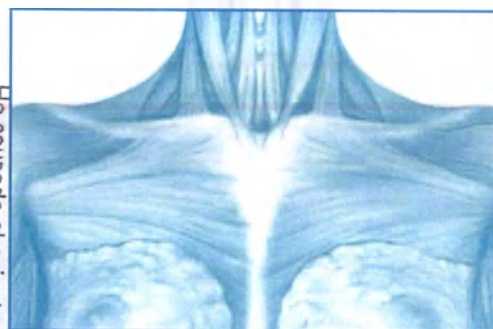
Conductive Cream Care Cream или
Conductive Cream Face Cream

RECOMMENDATIONS

It is recommended to carry out this procedure after targeted therapy or as concomitant therapy to improve the end result. Effectively combined with exercise and an appropriate diet.

COMMENTS

The proposed protocols are general and should be adapted to the individual requirements or circumstances of the patient, taking into account age, hormone imbalance, concomitant conditions, lifestyle, quality of life and history of surgical intervention. All of these factors can affect the end result. For external use only.


APPLICATION TECHNOLOGY

1. **Position:** hold the electrode evenly on the surface of the body; do not tilt.
2. **Pressure:** maintain sufficient pressure to ensure adequate contact between the skin and the surface of the electrode. The pressure should be stable and not cause discomfort to the patient.
3. **Speed:** move the electrode in a circular motion, each time slowly increasing the diameter of the circle. Movements should be slow (one circle per second) and smooth.
4. **Температура:** ensure that the temperature is suitable for the patient and is comfortable.
5. **Гигиена:** To optimize the transfer of energy to the patient when using the product, it is recommended that the user hold onto the electrode handle during the procedure. If it is necessary to use manual therapy methods, INDIBA representatives recommend the use of gloves (not included). This undoubtedly optimizes the transfer of energy to the patient and ensures the best result of the procedure.

CONTRAINDICATIONS

The procedure should not be carried out in the presence of a pacemaker or other electronic implants, in pregnant women, in case of skin wounds (open wounds or recent burns), thrombophlebitis and lack of sensitivity. Not contraindicated in the presence of implants

Patient sensations based on IAS value (IAS = INDIBA analog scale)

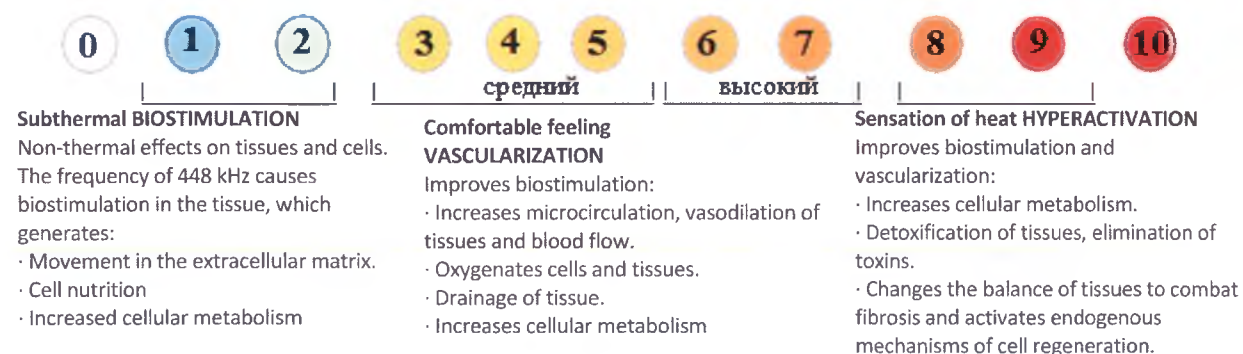CHEST (30 min)

Treatment with the aim of lifting, improving, tightening the neckline and reducing fibrosis.

POSITION

Lying on your back. Return flat electrode between the blades.

DESCRIPTION OF STAGES OF THE PROCEDURE

CAP 8 min IAS* 7

RES 5 min IAS* 7

EACH CHEST

CAP 2 min IAS* 2

APPLICATION DATA

Perform circular movements in the direction of the muscles and muscle cords. At the final stage (resistive mode) apply the methods of "tissue remodeling".

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Maintenance therapy: once a month

NOTES

If the use of cosmetics is necessary, they should be applied after the procedure

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Conductive Cream Care Cream

RECOMMENDATIONS

It is recommended to carry out this procedure after targeted therapy or as concomitant therapy to improve the end result. Effectively combined with exercise and an appropriate diet.

COMMENTS

The proposed protocols are general and should be adapted to the individual requirements or circumstances of the patient, taking into account age, hormone imbalance, concomitant conditions, lifestyle, quality of life and history of surgical intervention. All of these factors can affect the end result. For external use only.


APPLICATION TECHNOLOGY

- 1. Position:** hold the electrode evenly on the surface of the body; do not tilt.
- 2. Pressure:** maintain sufficient pressure to ensure adequate contact between the skin and the surface of the electrode. The pressure should be stable and not cause discomfort to the patient.
- 3. Speed:** move the electrode in a circular motion, each time slowly increasing the diameter of the circle. Movements should be slow (one circle per second) and smooth.
- 4. Температура:** ensure that the temperature is suitable for the patient and is comfortable.
- 5. Гигиена:** To optimize the transfer of energy to the patient when using the product, it is recommended that the user hold onto the electrode handle during the procedure. If it is necessary to use manual therapy methods, INDIBA representatives recommend the use of gloves (not included). This undoubtedly optimizes the transfer of energy to the patient and ensures the best result of the procedure.

CONTRAINDICATIONS

The procedure should not be carried out in the presence of a pacemaker or other electronic implants, in pregnant women, in case of skin wounds (open wounds or recent burns), thrombophlebitis and lack of sensitivity. Not contraindicated in the presence of implants

Patient sensations based on IAS value (IAS = INDIBA analog scale)

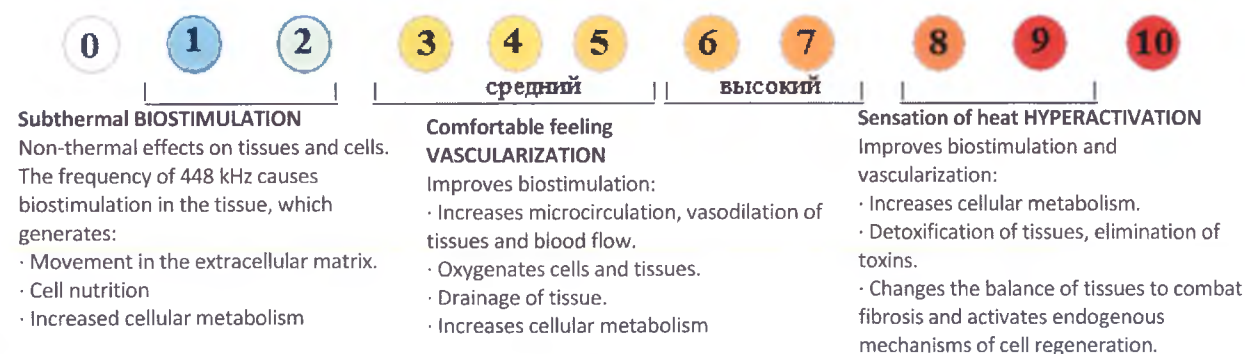ARMS (30 min)

INDIBA
DEEP CARE
MEDICAL FROENIC SYSTEM

Procedure for lightening the skin, reducing fat, shaping the contour of the hands and drainage.

POSITION

Lying on your back. Return flat electrode between the blades.

DESCRIPTION OF STAGES OF THE PROCEDURE

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APPLICATION DATA

Circular movements with an emphasis on areas with a lot of fat

SESSIONS

Two times a week

Total: 6 sessions

Maintenance therapy: once a month

NOTES

If the use of cosmetics is necessary, they should be applied after the procedure.

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Conductive Cream Care Cream

RECOMMENDATIONS

It is recommended to carry out this procedure after targeted therapy or as concomitant therapy to improve the end result. Effectively combined with exercise and an appropriate diet.

COMMENTS

The proposed protocols are general and should be adapted to the individual requirements or circumstances of the patient, taking into account age, hormone imbalance, concomitant conditions, lifestyle, quality of life and history of surgical intervention. All of these factors can affect the end result. For external use only.


APPLICATION TECHNOLOGY

- 1. Position:** hold the electrode evenly on the surface of the body; do not tilt.
- 2. Pressure:** maintain sufficient pressure to ensure adequate contact between the skin and the surface of the electrode. The pressure should be stable and not cause discomfort to the patient.
- 3. Speed:** move the electrode in a circular motion, each time slowly increasing the diameter of the circle. Movements should be slow (one circle per second) and smooth.
- 4. Температура:** ensure that the temperature is suitable for the patient and is comfortable.
- 5. Гигиена:** To optimize the transfer of energy to the patient when using the product, it is recommended that the user hold onto the electrode handle during the procedure. If it is necessary to use manual therapy methods, INDIBA representatives recommend the use of gloves (not included). This undoubtedly optimizes the transfer of energy to the patient and ensures the best result of the procedure.

CONTRAINDICATIONS

The procedure should not be carried out in the presence of a pacemaker or other electronic implants, in pregnant women, in case of skin wounds (open wounds or recent burns), thrombophlebitis and lack of sensitivity. Not contraindicated in the presence of implants

Patient sensations based on IAS value (IAS = INDIBA analog scale)

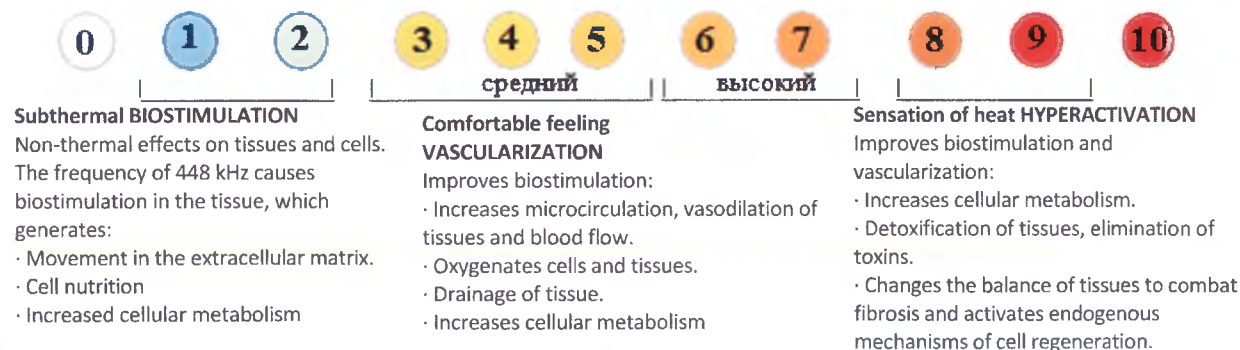ABDOMEN (30 min)

Procedure to reduce accumulated fat, strengthen and increase skin elasticity.

POSITION

Lying on your back. Return flat electrode between the blades.

DESCRIPTION OF STAGES OF THE PROCEDURE

CAP 16 min IAS* 8

RES 10 min IAS* 8

RES 4 min IAS* 2

APPLICATION DATA

Circular movements in areas of predominant accumulation of body fat. At the final stage (resistive mode), perform lymphatic drainage massage of the abdomen area.

SESSIONS

Two to three times a week


Total: 12 sessions

Maintenance therapy: once a month

NOTES

If the use of cosmetics is necessary, they should be applied after the procedure.

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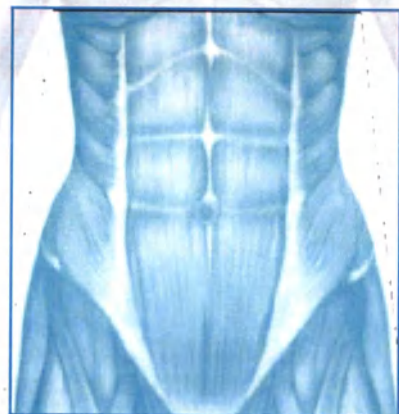RECOMMENDATIONS

It is recommended to carry out this procedure after targeted therapy or as concomitant therapy to improve the end result. Effectively combined with exercise and an appropriate diet.

COMMENTS

The proposed protocols are general and should be adapted to the individual requirements or circumstances of the patient, taking into account age, hormone imbalance, concomitant conditions, lifestyle, quality of life and history of surgical intervention. All of these factors can affect the end result. For external use only.


APPLICATION TECHNOLOGY

1. **Position:** hold the electrode evenly on the surface of the body; do not tilt.
2. **Pressure:** maintain sufficient pressure to ensure adequate contact between the skin and the surface of the electrode. The pressure should be stable and not cause discomfort to the patient.
3. **Speed:** move the electrode in a circular motion, each time slowly increasing the diameter of the circle. Movements should be slow (one circle per second) and smooth.
4. **Температура:** ensure that the temperature is suitable for the patient and is comfortable.
5. **Гигиена:** To optimize the transfer of energy to the patient when using the product, it is recommended that the user hold onto the electrode handle during the procedure. If it is necessary to use manual therapy methods, INDIBA representatives recommend the use of gloves (not included). This undoubtedly optimizes the transfer of energy to the patient and ensures the best result of the procedure.

CONTRAINDICATIONS

The procedure should not be carried out in the presence of a pacemaker or other electronic implants, in pregnant women, in case of skin wounds (open wounds or recent burns), thrombophlebitis and lack of sensitivity. Not contraindicated in the presence of implants

Patient sensations based on IAS value (IAS = INDIBA analog scale)

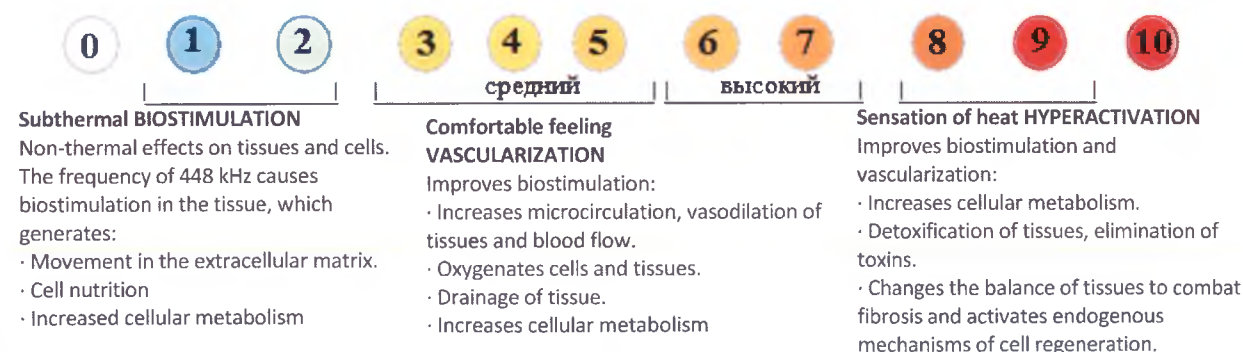BUTTOCKS (30 min)

A procedure to improve muscle and vascular tone, strengthen the skin, reduce localized fat, improve the contour, raise the buttocks, treat elastosis and sagging.

POSITION

Lying on the abdomen. Return flat electrode on the abdomen area

DESCRIPTION OF STAGES OF THE PROCEDURE

CAP 5 min IAS* 8 | EACH BUTTOCK

RES 8 min IAS* 8 | EACH BUTTOCK

RES 4 min IAS* 2 | SACRUM

APPLICATION DATA

Circular movements with an emphasis on areas of accumulation of body fat. At the final stage (resistive mode), go to the sacral region.

SESSIONS

Two to three times a week


Total: 12 sessions

Maintenance therapy: once a month

NOTES

If the use of cosmetics is necessary, they should be applied after the procedure.

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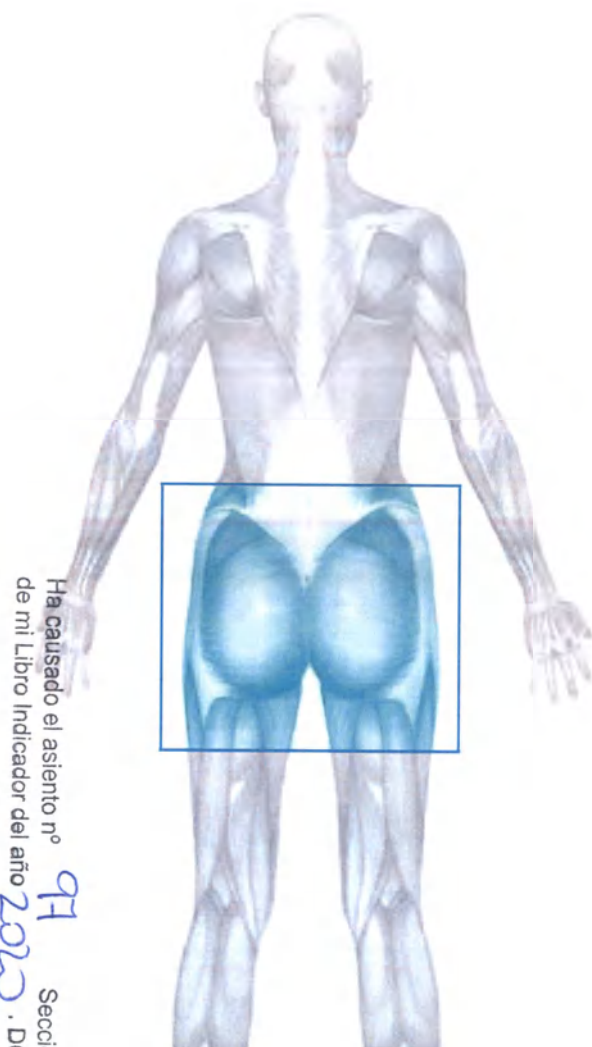RECOMMENDATIONS

It is recommended to carry out this procedure after targeted therapy or as concomitant therapy to improve the end result. Effectively combined with exercise and an appropriate diet.

COMMENTS

The proposed protocols are general and should be adapted to the individual requirements or circumstances of the patient, taking into account age, hormone imbalance, concomitant conditions, lifestyle, quality of life and history of surgical intervention. All of these factors can affect the end result. For external use only.


APPLICATION TECHNOLOGY

- 1. Position:** hold the electrode evenly on the surface of the body; do not tilt.
- 2. Pressure:** maintain sufficient pressure to ensure adequate contact between the skin and the surface of the electrode. The pressure should be stable and not cause discomfort to the patient.
- 3. Speed:** move the electrode in a circular motion, each time slowly increasing the diameter of the circle. Movements should be slow (one circle per second) and smooth.
- 4. Температура:** ensure that the temperature is suitable for the patient and is comfortable.
- 5. Гигиена:** To optimize the transfer of energy to the patient when using the product, it is recommended that the user hold onto the electrode handle during the procedure. If it is necessary to use manual therapy methods, INDIBA representatives recommend the use of gloves (not included). This undoubtedly optimizes the transfer of energy to the patient and ensures the best result of the procedure.

CONTRAINDICATIONS

The procedure should not be carried out in the presence of a pacemaker or other electronic implants, in pregnant women, in case of skin wounds (open wounds or recent burns), thrombophlebitis and lack of sensitivity. Not contraindicated in the presence of implants

Patient sensations based on IAS value (IAS = INDIBA analog scale)

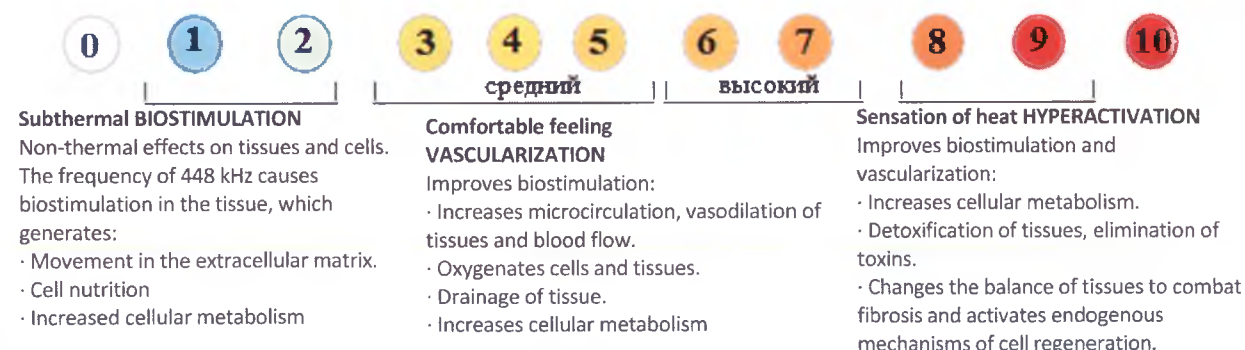This protocol, which is exclusively owned by INDIBA, S.A., is protected by intellectual property law and may only be used by the rightful owners of the relevant equipment. Therefore, it is forbidden to reproduce, modify, distribute or transmit it in whole or in part without the written consent of INDIBA, S.A. Все нарушения будут преследоваться по закону.

CELLULITE INDIBA[®] DEEP CARE MEDICAL PROIONIC SYSTEM

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acru

CELLULITE

The first part - lying on your back (30 min)

INDIBA
DEEP CARE
MEDICAL PROIONIC SYSTEM

Procedure to improve venous circulation, improve microcirculation, reduce swelling and increase skin elasticity

POSITION

Lying on the back: a return flat electrode on the lumbosacral region

DESCRIPTION OF STAGES OF THE PROCEDURE

CAP 2 min IAS* 6 | ABDOMEN
CAP 6 min IAS* 8
RES 3 min IAS* 8 | EACH LEG
CAP 2 min IAS* 2

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APPLICATION DATA

Circular movements with an emphasis on areas of accumulation of body fat. From the fourth session, change the time of the CAP and RES modes, taking into account the type of cellulite.

SESSIONS

Two times a week

Total: 15-20 sessions

NOTES

In case of edematous cellulitis, apply IAS 4 locally in the legs

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RECOMMENDATIONS

It is recommended to carry out this procedure after targeted therapy or as concomitant therapy to improve the end result. Effectively combined with exercise and an appropriate diet.

CONTRAINDICATIONS

The procedure should not be carried out in the presence of a pacemaker or other electronic implants, in pregnant women, in case of skin wounds (open wounds or recent burns), thrombophlebitis and lack of sensitivity. Not contraindicated in the presence of implants


COMMENTS

The proposed protocols are general and should be adapted to the individual requirements or circumstances of the patient, taking into account age, hormone imbalance, concomitant conditions, lifestyle, quality of life and history of surgical intervention. All of these factors can affect the end result. For external use only.

APPLICATION TECHNOLOGY

1. **Position:** hold the electrode evenly on the surface of the body; do not tilt.
2. **Pressure:** maintain sufficient pressure to ensure adequate contact between the skin and the surface of the electrode. The pressure should be stable and not cause discomfort to the patient.
3. **Speed:** move the electrode in a circular motion, each time slowly increasing the diameter of the circle. Movements should be slow (one circle per second) and smooth.
4. **Температура:** ensure that the temperature is suitable for the patient and is comfortable.
5. **Гигиена:** To optimize the transfer of energy to the patient when using the product, it is recommended that the user hold onto the electrode handle during the procedure. If it is necessary to use manual therapy methods, INDIBA representatives recommend the use of gloves (not included). This undoubtedly optimizes the transfer of energy to the patient and ensures the best result of the procedure.

Patient sensations based on IAS value (IAS = INDIBA analog scale)

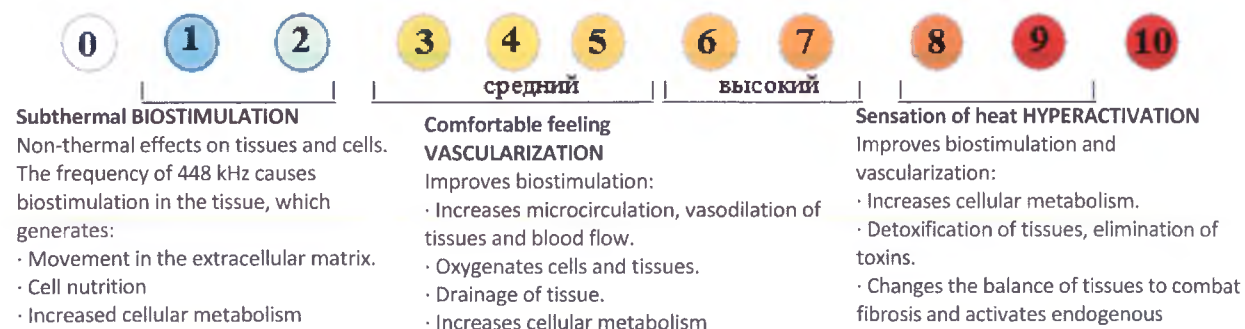
This protocol, which is exclusively owned by INDIBA, S.A., is protected by intellectual property law and may only be used by the rightful owners of the relevant equipment. Therefore, it is forbidden to reproduce, modify, distribute or transmit it in whole or in part without the written consent of INDIBA, S.A., issued in written form. All violations will be prosecuted by law.

CELLULITE

The second part - lying on the ABDOMEN (24 min)

INDIBA
DEEP CARE
MEDICAL PROIONIC SYSTEM

Procedure to improve venous circulation, improve microcirculation, reduce swelling and increase skin elasticity

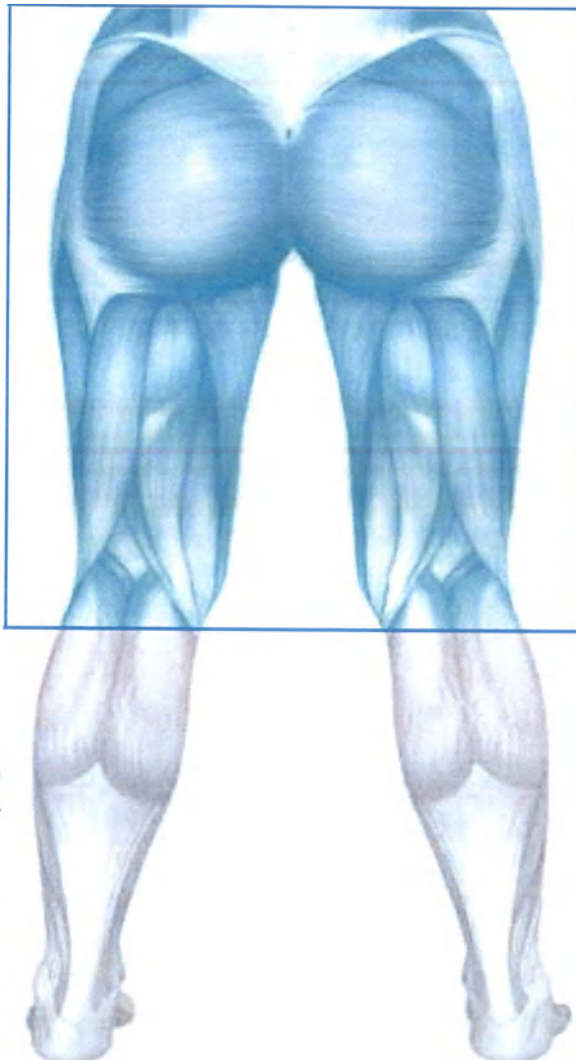POSITION

Lying on the ABDOMEN: return flat electrode on the ABDOMEN area

DESCRIPTION OF STAGES OF THE PROCEDURE

CAP 4 min IAS* 8		
CAP 4 min IAS* 8		EVERY LEG
RES 3 min IAS* 2		
CAP 2 min IAS* 4		SACRUM

APPLICATION DATA

Circular movements with an emphasis on areas of accumulation of body fat. From the fourth session, change the time of the CAP and RES modes, taking into account the CELLULITE type.

SESSIONS

Two times a week

Total: 15-20 sessions

NOTES

In case of edematous CELLULITE, apply IAS 4 locally in the leg area

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CELLULITE

The second part - lying on the ABDOMEN

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RECOMMENDATIONS

It is recommended to carry out this procedure after targeted therapy or as concomitant therapy to improve the end result. Effectively combined with exercise and an appropriate diet.

COMMENTS

The proposed protocols are general and should be adapted to the individual requirements or circumstances of the patient, taking into account age, hormone imbalance, concomitant conditions, lifestyle, quality of life and history of surgical intervention. All of these factors can affect the end result. For external use only.


APPLICATION TECHNOLOGY

1. **Position:** hold the electrode evenly on the surface of the body; do not tilt.
2. **Pressure:** maintain sufficient pressure to ensure adequate contact between the skin and the surface of the electrode. The pressure should be stable and not cause discomfort to the patient.
3. **Speed:** move the electrode in a circular motion, each time slowly increasing the diameter of the circle. Movements should be slow (one circle per second) and smooth.
4. **Температура:** ensure that the temperature is suitable for the patient and is comfortable.
5. **Гигиена:** To optimize the transfer of energy to the patient when using the product, it is recommended that the user hold onto the electrode handle during the procedure. If it is necessary to use manual therapy methods, INDIBA representatives recommend the use of gloves (not included). This undoubtedly optimizes the transfer of energy to the patient and ensures the best result of the procedure.

CONTRAINDICATIONS

The procedure should not be carried out in the presence of a pacemaker or other electronic implants, in pregnant women, in case of skin wounds (open wounds or recent burns), thrombophlebitis and lack of sensitivity. Not contraindicated in the presence of implants

Patient sensations based on IAS value (IAS = INDIBA analog scale)

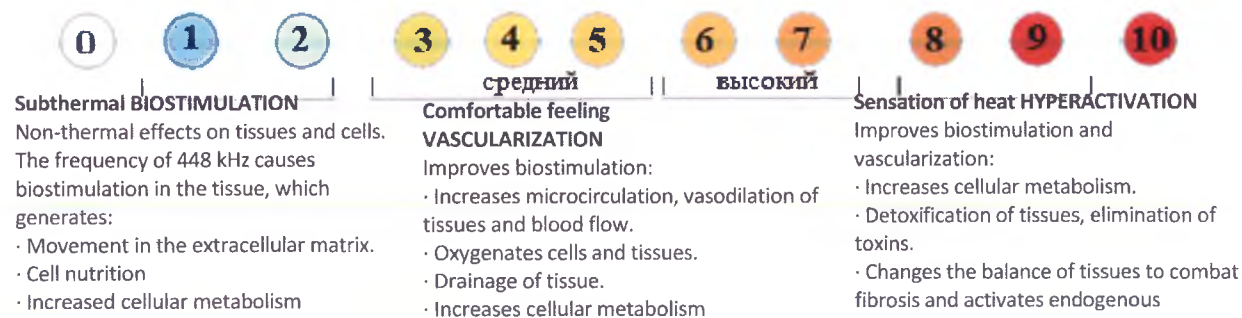This protocol, which is exclusively owned by INDIBA, S.A., is protected by intellectual property law and may only be used by the rightful owners of the relevant equipment. Therefore, it is forbidden to reproduce, modify, distribute or transmit it in whole or in part without the written consent of INDIBA, S.A., issued in written form. All violations will be prosecuted by law. *

OTHER APPLICATIONS

INDIBA
DEEP CARE
MEDICAL PROIONIC SYSTEM

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MASSAGE NOZZLE

INDIBA
DEEP CARE
MEDICAL PROIONIC SYSTEM

The massage nozzle is an optional product for the 65 mm diameter RES electrode, which adds a mechanical effect to electrical and thermal effects. Due to this, the impulses penetrate more deeply during the procedure.

COMPATIBLE ELECTRODE

RES electrode with a diameter of 65 mm

OBJECTIVE OF THE PROCEDURE

- Improved drainage
- Reduced body fat
- Softening fibrous tissue
- Deeper massage without pain

APPLY FOR PROTOCOLS

- Before surgery (body and limbs)
- After surgery (post liposuction, part 2, prevention and fibrous tissue)
- Skin tightening procedures (ABDOMEN and limbs)
- Treatment of CELLULITE
- Anti-aging procedure for the back

APPLICATION DATA

Circular movements with an emphasis on areas of accumulation of fatty deposits and on fibrous tissues.

OTHER COMMENTS

Should be used after capacitive mode to prepare the fabric.

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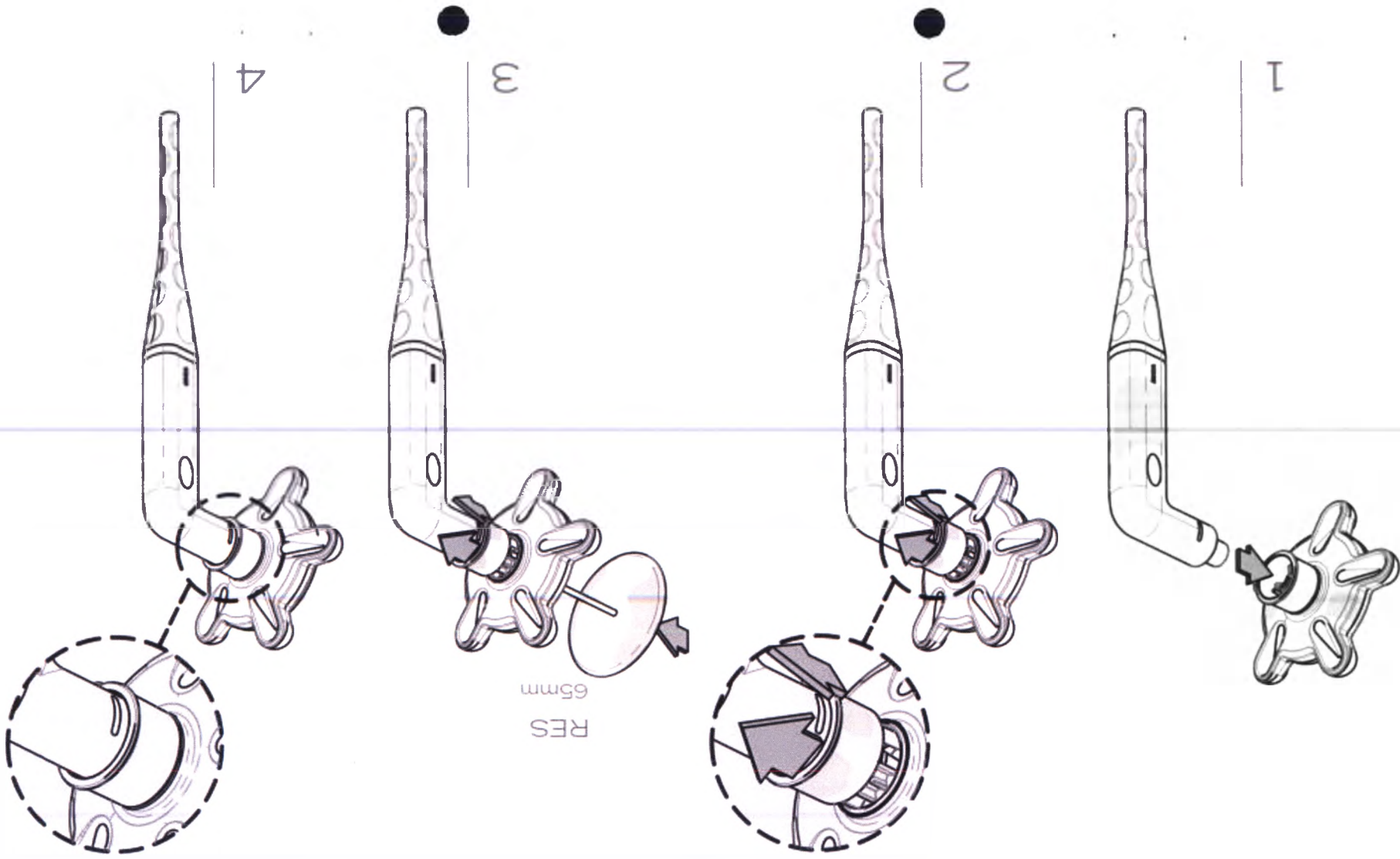
APPLICATION DATA

Linear, circular movements and drainage lines

SESSIONS


Two times a week

Total: 12 sessions

NOTES

It is recommended to do exercises for posture and muscle strengthening.

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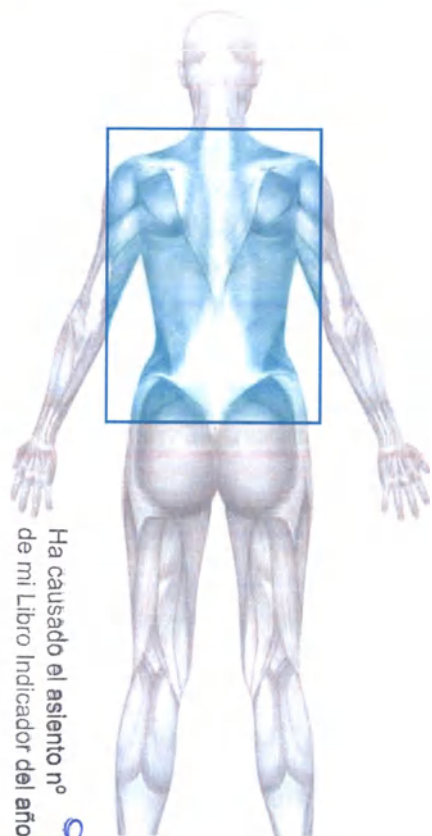BACK ANTI-AGING PROCEDURE

INDIBA
DEEP CARE
MEDICAL PROIONIC SYSTEM

USED PRODUCTS

Conductive Cream Care Cream

COMMENTS

The proposed protocols are general and should be adapted to the individual requirements or circumstances of the patient, taking into account age, hormone imbalance, concomitant conditions, lifestyle, quality of life and history of surgical intervention. All of these factors can affect the end result. For external use only.

APPLICATION TECHNOLOGY

1. **Position:** hold the electrode evenly on the surface of the body; do not tilt.
2. **Pressure:** maintain sufficient pressure to ensure adequate contact between the skin and the surface of the electrode. The pressure should be stable and not cause discomfort to the patient.
3. **Speed:** move the electrode in a circular motion, each time slowly increasing the diameter of the circle. Movements should be slow (one circle per second) and smooth.
4. **Температура:** ensure that the temperature is suitable for the patient and is comfortable.
5. **Гигиена:** To optimize the transfer of energy to the patient when using the product, it is recommended that the user hold onto the electrode handle during the procedure. If it is necessary to use manual therapy methods, INDIBA representatives recommend the use of gloves (not included). This undoubtedly optimizes the transfer of energy to the patient and ensures the best result of the procedure.


RECOMMENDATIONS

It is recommended to carry out this procedure after targeted therapy or as concomitant therapy to improve the end result. Effectively combined with exercise and an appropriate diet.

CONTRAINDICATIONS

The procedure should not be carried out in the presence of a pacemaker or other electronic implants, in pregnant women, in case of skin wounds (open wounds or recent burns), thrombophlebitis and lack of sensitivity. Not contraindicated in the presence of implants

Patient sensations based on IAS value (IAS = INDIBA analog scale)

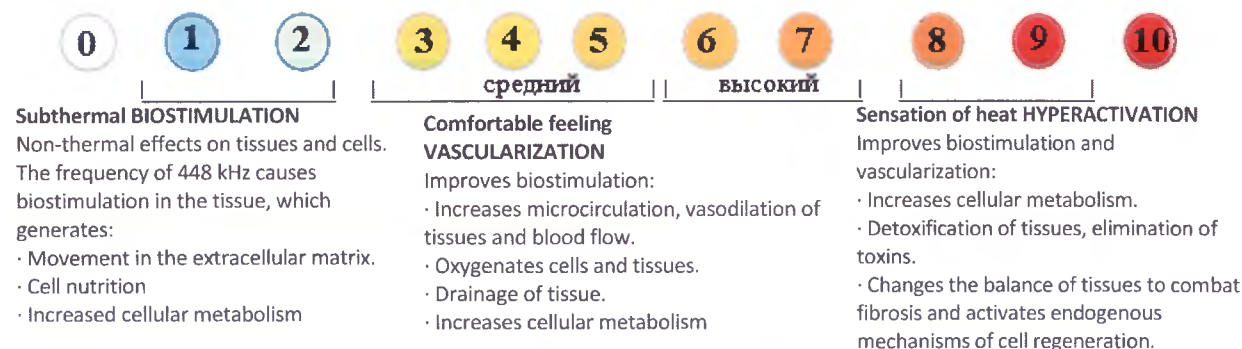This protocol, which is exclusively owned by INDIBA, S.A., is protected by intellectual property law and may only be used by the rightful owners of the relevant equipment. Therefore, it is forbidden to reproduce, modify, distribute or transmit it in whole or in part without the prior written consent of INDIBA, S.A. All violations will be prosecuted by law.

PREVENTION OF DEEP VEIN THROMBOSIS (40 min)

INDIBA
DEEP CARE
MEDICAL PROIONIC SYSTEM

Improvement of venous circulation, microcirculation, pain and drainage control.

POSITION

Lying on the back, return flat electrode on the lumbosacral region

DESCRIPTION OF STAGES OF THE PROCEDURE

CAP 7 min IAS * 5 (lying on your back and ABDOMEN)

EACH LEG

RES 2.5 min IAS * 5 (lying on the back ABDOMEN, fixed on the sole)

RES 1 min IAS * 3 (fixed on the sole of the foot)

APPLICATION DATA

Circular and linear movements in the cranial direction

SESSIONS

Two times a week

Total: **12 sessions**

NOTES

It is recommended to use cosmetics for vascular diseases. Especially indicated for people with standing work

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97

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PREVENTION OF DEEP VEIN THROMBOSIS

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RECOMMENDATIONS

It is recommended to carry out this procedure after targeted therapy or as concomitant therapy to improve the end result. Effectively combined with exercise and an appropriate diet.

COMMENTS

The proposed protocols are general and should be adapted to the individual requirements or circumstances of the patient, taking into account age, hormone imbalance, concomitant conditions, lifestyle, quality of life and history of surgical intervention. All of these factors can affect the end result. For external use only.


APPLICATION TECHNOLOGY

- 1. Position:** hold the electrode evenly on the surface of the body; do not tilt.
- 2. Pressure:** maintain sufficient pressure to ensure adequate contact between the skin and the surface of the electrode. The pressure should be stable and not cause discomfort to the patient.
- 3. Speed:** move the electrode in a circular motion, each time slowly increasing the diameter of the circle. Movements should be slow (one circle per second) and smooth.
- 4. Температура:** ensure that the temperature is suitable for the patient and is comfortable.
- 5. Гигиена:** To optimize the transfer of energy to the patient when using the product, it is recommended that the user hold onto the electrode handle during the procedure. If it is necessary to use manual therapy methods, INDIBA representatives recommend the use of gloves (not included). This undoubtedly optimizes the transfer of energy to the patient and ensures the best result of the procedure.

CONTRAINDICATIONS

The procedure should not be carried out in the presence of a pacemaker or other electronic implants, in pregnant women, in case of skin wounds (open wounds or recent burns), thrombophlebitis and lack of sensitivity. Not contraindicated in the presence of implants

Patient sensations based on IAS value (IAS = INDIBA analog scale)

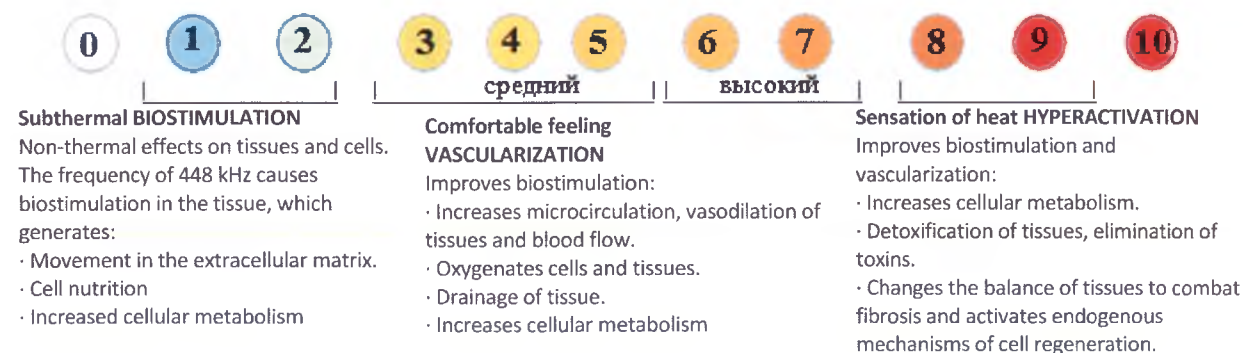TREATMENT OF STRETCHES (30 min)

Improving vascularization, pain control and drainage.

POSITION

The return flat electrode in the interscapular or lumbar region in the supine position (depending on the area of impact)

CAP 15 min IAS* 7

RES 15 min IAS* 7

DESCRIPTION OF STAGES OF THE PROCEDURE

APPLICATION DATA

Local application in the area requiring treatment

SESSIONS

Once a week

Total: 4 sessions

NOTES

Repeat the procedure after a month and a half using the following steps:

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RECOMMENDATIONS

It is recommended to carry out this procedure after targeted therapy or as concomitant therapy to improve the end result. Effectively combined with exercise and an appropriate diet.

COMMENTS

The proposed protocols are general and should be adapted to the individual requirements or circumstances of the patient, taking into account age, hormone imbalance, concomitant conditions, lifestyle, quality of life and history of surgical intervention. All of these factors can affect the end result. For external use only.


APPLICATION TECHNOLOGY

- 1. Position:** hold the electrode evenly on the surface of the body; do not tilt.
- 2. Pressure:** maintain sufficient pressure to ensure adequate contact between the skin and the surface of the electrode. The pressure should be stable and not cause discomfort to the patient.
- 3. Speed:** move the electrode in a circular motion, each time slowly increasing the diameter of the circle. Movements should be slow (one circle per second) and smooth.
- 4. Температура:** ensure that the temperature is suitable for the patient and is comfortable.
- 5. Гигиена:** To optimize the transfer of energy to the patient when using the product, it is recommended that the user hold onto the electrode handle during the procedure. If it is necessary to use manual therapy methods, INDIBA representatives recommend the use of gloves (not included). This undoubtedly optimizes the transfer of energy to the patient and ensures the best result of the procedure.


CONTRAINDICATIONS

The procedure should not be carried out in the presence of a pacemaker or other electronic implants, in pregnant women, in case of skin wounds (open wounds or recent burns), thrombophlebitis and lack of sensitivity. Not contraindicated in the presence of implants

Patient sensations based on IAS value (IAS = INDIBA analog scale)

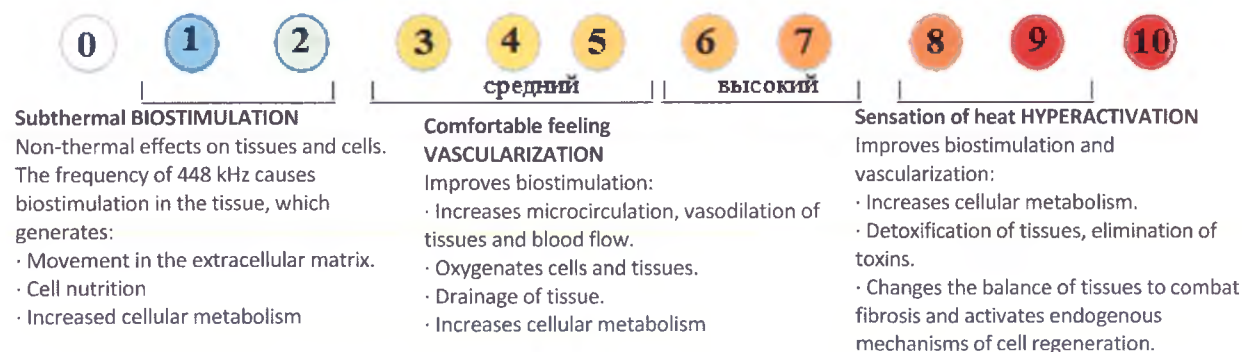This protocol, which is exclusively owned by INDIBA, S.A., is protected by intellectual property law and may only be used by the rightful owners of the relevant equipment. Therefore, it is forbidden to reproduce, modify, distribute or transmit it in whole or in part without the written consent of INDIBA, S.A., issued in written form. All violations will be pursued by law.

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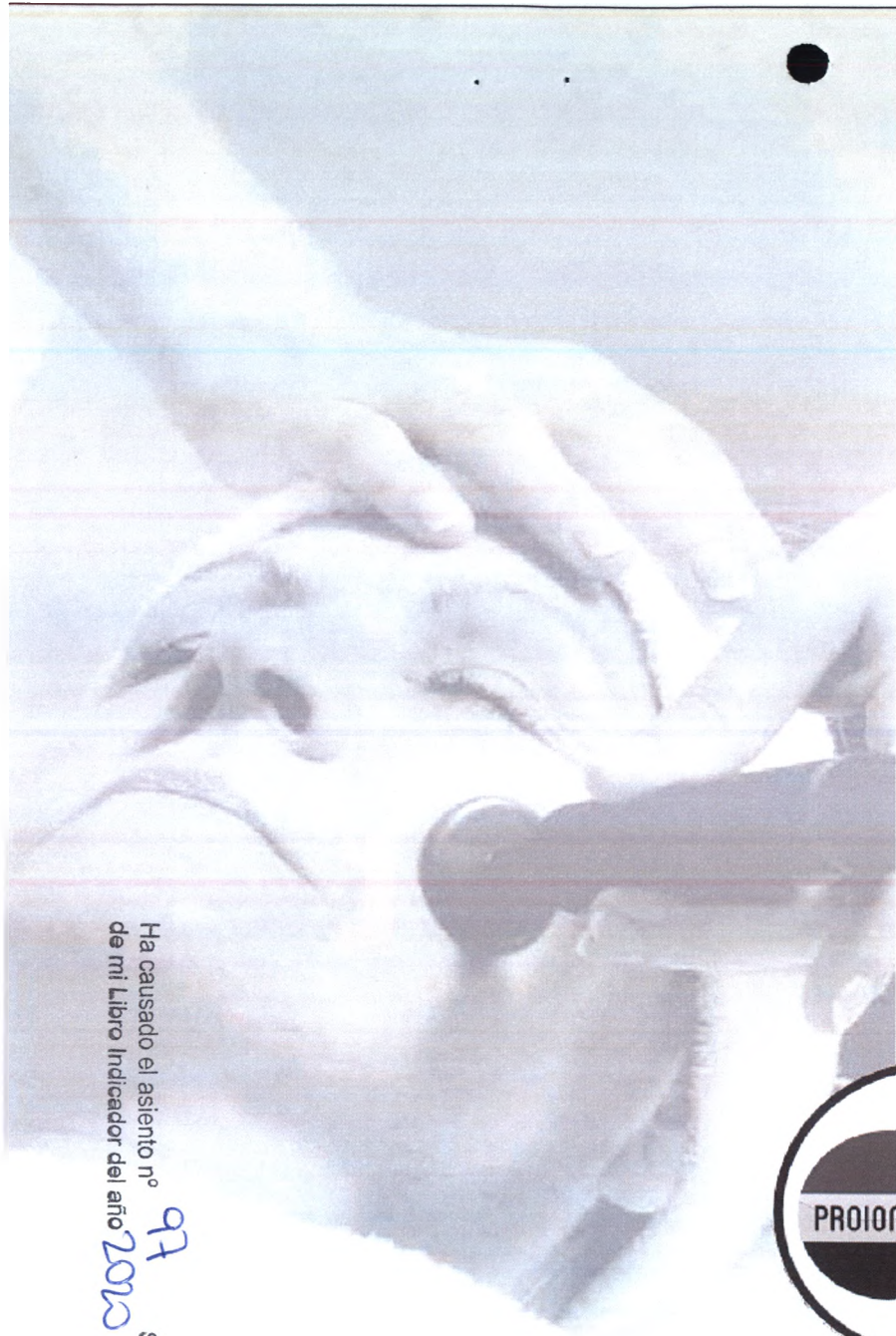COMBINED PROTOCOLS

INDIBA
DEEP CARE
MEDICAL PROIONIC SYSTEM

COMBINED PROTOCOLS

BEFORE THE FILING PROCEDURE PDO (polydioxanone), PCL (polycaprolactone), and PLL (polylactic) (FACE AND NECK 30 min)

INDIBA
DEEP CARE
MEDICAL PROIONIC SYSTEM

Preparation of tissue for treatment, increased blood flow and oxygenation, contributing to a higher therapeutic effect due to increased cellular metabolism when installing an implant from threads. In addition, the doctor will feel less resistance during the placement of threads, and the risk of hematomas will decrease

POSITION

Lying on your back. Return flat electrode on the interscapular region (ensure good contact with the tissue)

DESCRIPTION OF STAGES OF THE PROCEDURE

CAP 20 - 30 min IAS* 6 - 8

APPLICATION DATA

ELECTRODE CAP: 30 - 40 mm

Move along the drainage line and towards the musculature

SESSIONS

3 sessions. 30 - 15 and 1 day before implantation of threads

Total: 3 sessions

BEGINNING OF THE PROCEDURE

30 days before implantation of the threads

NOTES

The procedure can be carried out once or as an adjunct therapy to any other treatment determined by the doctor, for example, light-emitting diodes (LEDs)

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RECOMMENDATIONS

The installation of threads is recommended to be carried out 24 hours after the procedure in order to eliminate the thermal effect.

Check for past hypersensitivity to face cream ingredients Face Cream

COMMENTS

The proposed protocols are general and should be adapted to the individual requirements or circumstances of the patient, taking into account age, hormone imbalance, concomitant conditions, lifestyle, quality of life and history of surgical intervention. All of these factors can affect the end result. For external use only.


APPLICATION TECHNOLOGY

1. **Position:** hold the electrode evenly on the surface of the body; do not tilt.
2. **Pressure:** maintain sufficient pressure to ensure adequate contact between the skin and the surface of the electrode. The pressure should be stable and not cause discomfort to the patient.
3. **Speed:** move the electrode in a circular motion, each time slowly increasing the diameter of the circle. Movements should be slow (one circle per second) and smooth.
4. **Температура:** ensure that the temperature is suitable for the patient and is comfortable.
5. **Гигиена:** To optimize the transfer of energy to the patient when using the product, it is recommended that the user hold onto the electrode handle during the procedure. If it is necessary to use manual therapy methods, INDIBA representatives recommend the use of gloves (not included). This undoubtedly optimizes the transfer of energy to the patient and ensures the best result of the procedure.

CONTRAINDICATIONS

The procedure should not be carried out in the presence of a pacemaker or other electronic implants, in pregnant women, in case of skin wounds (open wounds or recent burns), thrombophlebitis and lack of sensitivity. Not contraindicated in the presence of implants

Patient sensations based on IAS value (IAS = INDIBA analog scale)

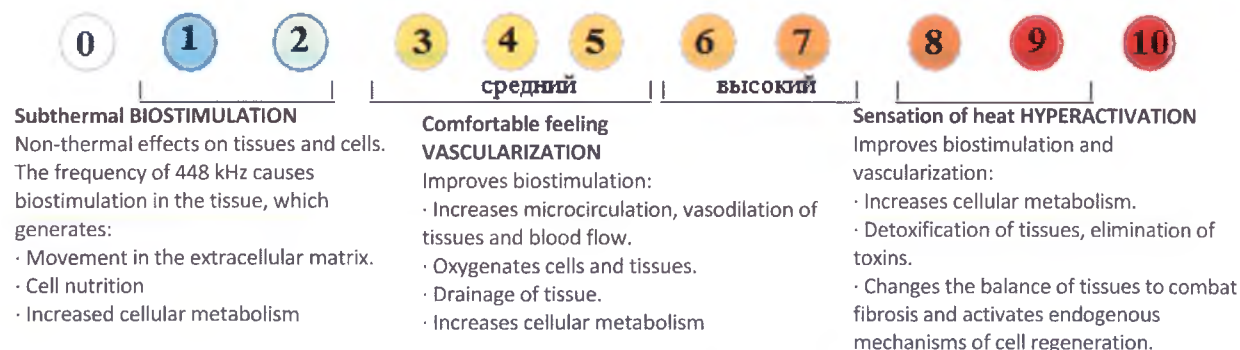This protocol, which is exclusively owned by INDIBA, S.A., is protected by intellectual property law and may only be used by the rightful owners of the relevant equipment. Therefore, it is forbidden to reproduce, modify, distribute or transmit it in whole or in part without the written consent of INDIBA, S.A., issued in written form. All violations will be prosecuted by law.

AFTER THE FILING PROCEDURE PDO (polydioxanone), PCL (polycaprolactone) and PLL (polylactic) (FACE AND NECK 20-30 min)

INDIBA
DEEP CARE
MEDICAL PROIONIC SYSTEM

Strengthening neocollagenesis caused by the placement of the thread itself, and contributing to the formation of physiological reactions. Improving the results.

POSITION

Lying on your back. Return flat electrode on the interscapular region (ensure good contact with the tissue)

DESCRIPTION OF STAGES OF THE PROCEDURE

CAP 10 min **IAS*** 6 – 8

RES 10-20 min **IAS*** 6 – 8

ELECTRODE CAP: 30 - 40 mm, **RES** 35 mm

APPLICATION DATA

Move along the drainage line and towards the musculature

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NOTES

If necessary, a subthermal procedure (IAS 1 with signal modulation) can be performed immediately after placement of the threads to limit the side effects of the operation at the subcutaneous level, but without reducing the effectiveness of the treatment. This protocol will consist of 1 to 3 sessions with an interval of 24 hours between sessions in accordance with the needs of the patient.

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BEFORE THE FILING PROCEDURE PDO (polydioxanone), PCL (polycaprolactone), and PLL (polylactic)

INDIBA
DEEP CARE
MEDICAL FREQUENT SYSTEM

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COMMENTS

The proposed protocols are general and should be adapted to the individual requirements or circumstances of the patient, taking into account age, hormone imbalance, concomitant conditions, lifestyle, quality of life and history of surgical intervention. All of these factors can affect the end result. For external use only.


APPLICATION TECHNOLOGY

- 1. Position:** hold the electrode evenly on the surface of the body; do not tilt.
- 2. Pressure:** maintain sufficient pressure to ensure adequate contact between the skin and the surface of the electrode. The pressure should be stable and not cause discomfort to the patient.
- 3. Speed:** move the electrode in a circular motion, each time slowly increasing the diameter of the circle. Movements should be slow (one circle per second) and smooth.
- 4. Температура:** ensure that the temperature is suitable for the patient and is comfortable.
- 5. Гигиена:** To optimize the transfer of energy to the patient when using the product, it is recommended that the user hold onto the electrode handle during the procedure. If it is necessary to use manual therapy methods, INDIBA representatives recommend the use of gloves (not included). This undoubtedly optimizes the transfer of energy to the patient and ensures the best result of the procedure.

CONTRAINDICATIONS

The procedure should not be carried out in the presence of a pacemaker or other electronic implants, in pregnant women, in case of skin wounds (open wounds or recent burns), thrombophlebitis and lack of sensitivity. Not contraindicated in the presence of implants

Patient sensations based on IAS value (IAS = INDIBA analog scale)

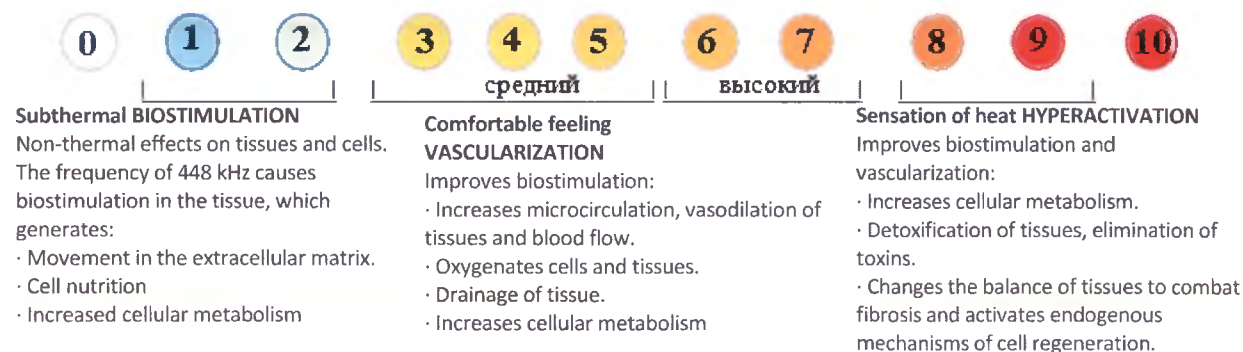This protocol, which is exclusively owned by INDIBA, S.A., is protected by intellectual property law and may only be used by the rightful owners of the relevant equipment. Therefore, it is forbidden to reproduce, modify, distribute or transmit it in whole or in part without the written consent of INDIBA, S.A., issued in written form. All violations will be prosecuted by law.

BEFORE INTRODUCING HYALURONIC ACID

(FACE 15 min)

INDIBA
DEEP CARE
MEDICAL PROIONIC SYSTEM

Immediate densification of the tissue, improved non-dispersion of the product caused by sagging skin of the patient, and, consequently, an increase in the volume of the treated area with a smaller amount of product.

POSITION

Lying on your back. Return flat electrode on the interscapular region (ensure good contact with the tissue)

DESCRIPTION OF STAGES OF THE PROCEDURE

CAP 10 min IAS* 6 – 8

RES 5 min IAS* 6 – 8

ELECTRODE CAP: 30 - 40 mm, RES 35 mm

APPLICATION DATA

The procedure is carried out throughout the face, focusing on the areas where it is planned to introduce hyaluronic acid

SESSIONS

Two or three times a week

Total: 2 or 4 sessions

BEGINNING OF THE PROCEDURE

1-2 weeks before administration

NOTES

If necessary, chemical peeling can be performed as an additional treatment 24/48 hours before the procedure.

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RECOMMENDATIONS

Drink at least 1.5 liters of water per day while maintaining a healthy lifestyle.

Consider a possible decrease in the therapeutic effectiveness of this procedure in patients with collagen diseases or patients taking high doses of corticosteroids.

Check for past hypersensitivity to face cream ingredients Face Cream

COMMENTS

The proposed protocols are general and should be adapted to the individual requirements or circumstances of the patient, taking into account age, hormone imbalance, concomitant conditions, lifestyle, quality of life and history of surgical intervention. All of these factors can affect the end result. For external use only.


APPLICATION TECHNOLOGY

- 1. Position:** hold the electrode evenly on the surface of the body; do not tilt.
- 2. Pressure:** maintain sufficient pressure to ensure adequate contact between the skin and the surface of the electrode. The pressure should be stable and not cause discomfort to the patient.
- 3. Speed:** move the electrode in a circular motion, each time slowly increasing the diameter of the circle. Movements should be slow (one circle per second) and smooth.
- 4. Температура:** ensure that the temperature is suitable for the patient and is comfortable.
- 5. Гигиена:** To optimize the transfer of energy to the patient when using the product, it is recommended that the user hold onto the electrode handle during the procedure. If it is necessary to use manual therapy methods, INDIBA representatives recommend the use of gloves (not included). This undoubtedly optimizes the transfer of energy to the patient and ensures the best result of the procedure.

CONTRAINDICATIONS

The procedure should not be carried out in the presence of a pacemaker or other electronic implants, in pregnant women, in case of skin wounds (open wounds or recent burns), thrombophlebitis and lack of sensitivity. Not contraindicated in the presence of implants

Patient sensations based on IAS value (IAS = INDIBA analog scale)

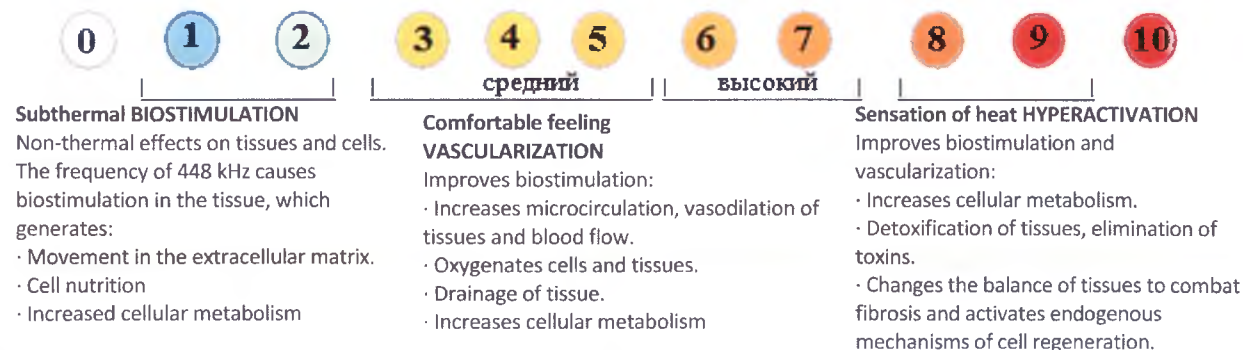This protocol, which is exclusively owned by INDIBA, S.A., is protected by intellectual property law and may only be used by the rightful owners of the relevant equipment. Therefore, it is forbidden to reproduce, modify, distribute or transmit it in whole or in part without the written consent of INDIBA, S.A. Все нарушения будут преследоваться по закону.

AFTER INTRODUCTION OF HYALURONIC ACID (FACE 15 min)

INDIBA
DEEP CARE
MEDICAL PROIONIC SYSTEM

Strengthening the biological effects of hyaluronic acid on the treated tissue (neocollagenesis, neoVASCULARIZATION, etc.). In addition, it is also possible to reduce the side effects of treatment (bruises, immobilization, edema ...).

POSITION

Lying on your back. Return flat electrode on the interscapular region (ensure good contact with the tissue)

DESCRIPTION OF STAGES OF THE PROCEDURE

CAP 10 min IAS* 3 – 5

RES 5 min IAS* 3 – 5

ELECTRODE CAP: 30 - 40 mm, RES 35 mm

APPLICATION DATA

Performing movements along the drainage line and towards the muscles

SESSIONS

One-Two times a week

Total: 2 or 4 sessions

BEGINNING OF THE PROCEDURE

15 days after administration

NOTES

In case of bruising, it is recommended to immediately begin treatment, leaving at least 8 hours between sessions, depending on the needs of the patient. Perform circular motions on bruises with the help of ELECTRODE-CAP 5 - 7 minutes with IAS 1 - 2 with signal modulation.

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Sección 2ª

[Handwritten signature]

AFTER INTRODUCTION OF HYALURONIC ACID

INDIBA
DEEP CARE
MEDICAL FREQUENT SYSTEM

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RECOMMENDATIONS

15 days after the end of mesotherapy, you can carry out the procedure as an adjuvant (IAS 3 - 5) against aging.

Check for past hypersensitivity to face cream ingredients Face Cream

COMMENTS

The proposed protocols are general and should be adapted to the individual requirements or circumstances of the patient, taking into account age, hormone imbalance, concomitant conditions, lifestyle, quality of life and history of surgical intervention. All of these factors can affect the end result. For external use only.


APPLICATION TECHNOLOGY

- 1. Position:** hold the electrode evenly on the surface of the body; do not tilt.
- 2. Pressure:** maintain sufficient pressure to ensure adequate contact between the skin and the surface of the electrode. The pressure should be stable and not cause discomfort to the patient.
- 3. Speed:** move the electrode in a circular motion, each time slowly increasing the diameter of the circle. Movements should be slow (one circle per second) and smooth.
- 4. Температура:** ensure that the temperature is suitable for the patient and is comfortable.
- 5. Гигиена:** To optimize the transfer of energy to the patient when using the product, it is recommended that the user hold onto the electrode handle during the procedure. If it is necessary to use manual therapy methods, INDIBA representatives recommend the use of gloves (not included). This undoubtedly optimizes the transfer of energy to the patient and ensures the best result of the procedure.

CONTRAINDICATIONS

The procedure should not be carried out in the presence of a pacemaker or other electronic implants, in pregnant women, in case of skin wounds (open wounds or recent burns), thrombophlebitis and lack of sensitivity. Not contraindicated in the presence of implants

Patient sensations based on IAS value (IAS = INDIBA analog scale)

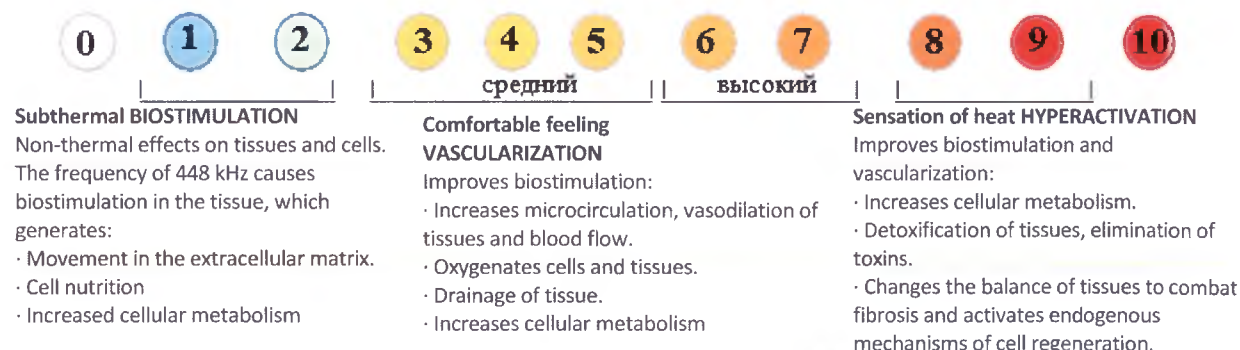
This protocol, which is exclusively owned by INDIBA, S.A., is protected by intellectual property law and may only be used by the rightful owners of the relevant equipment. Therefore, it is forbidden to reproduce, modify, distribute or transmit it in whole or in part without the written consent of INDIBA, S.A. Все нарушения будут преследоваться по закону.

BEFORE THE INTRODUCTION OF CALCIUM HYDROXYAPATITE, POLYCAPROLACTON OR POLYLACTIC ACID (FACE 20 min)

INDIBA
DEEP CARE
MEDICAL PROIONIC SYSTEM

Ensuring optimal tissue condition before infiltration, as a result of immediate tissue densification and increased cellular metabolism.

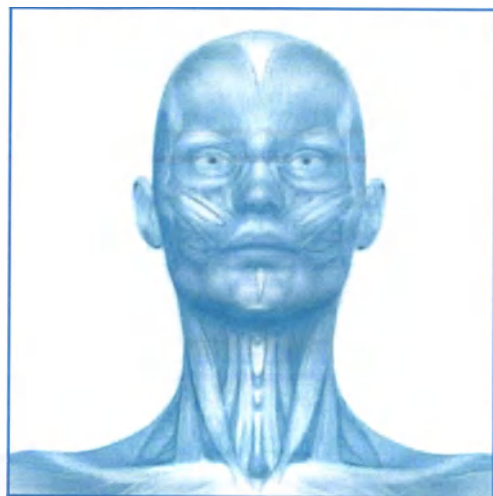POSITION

Lying on your back. Return flat electrode on the interscapular region (ensure good contact with the tissue)

DESCRIPTION OF STAGES OF THE PROCEDURE

CAP 15 min IAS* 6 – 8

RES 5 min IAS* 6 – 8

ELECTRODE CAP: 30 - 40 mm, RES 35 mm

APPLICATION DATA

The procedure is carried out throughout the face, focusing on the areas in which it is planned to introduce these substances

SESSIONS

Two times a week

Total: 4 sessions

BEGINNING OF THE PROCEDURE

15 days before administration

NOTES

The procedure can be carried out once or as an adjunct therapy to any other treatment determined by the doctor, for example, light-emitting diodes (LEDs)

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97

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BEFORE THE INTRODUCTION OF CALCIUM HYDROXYAPATITE, POLYCAPROLACTON OR POLYLACTIC ACID

INDIBA
DEEP CARE
MEDICAL FREQUENT SYSTEM

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RECOMMENDATIONS

The introduction of fillers or other specified substances is recommended 24 hours after the procedure to eliminate the thermal effect.

Check for past hypersensitivity to face cream ingredients Face Cream

COMMENTS

The proposed protocols are general and should be adapted to the individual requirements or circumstances of the patient, taking into account age, hormone imbalance, concomitant conditions, lifestyle, quality of life and history of surgical intervention. All of these factors can affect the end result. For external use only.


APPLICATION TECHNOLOGY

1. **Position:** hold the electrode evenly on the surface of the body; do not tilt.
2. **Pressure:** maintain sufficient pressure to ensure adequate contact between the skin and the surface of the electrode. The pressure should be stable and not cause discomfort to the patient.
3. **Speed:** move the electrode in a circular motion, each time slowly increasing the diameter of the circle. Movements should be slow (one circle per second) and smooth.
4. **Температура:** ensure that the temperature is suitable for the patient and is comfortable.
5. **Гигиена:** To optimize the transfer of energy to the patient when using the product, it is recommended that the user hold onto the electrode handle during the procedure. If it is necessary to use manual therapy methods, INDIBA representatives recommend the use of gloves (not included). This undoubtedly optimizes the transfer of energy to the patient and ensures the best result of the procedure.

CONTRAINDICATIONS

The procedure should not be carried out in the presence of a pacemaker or other electronic implants, in pregnant women, in case of skin wounds (open wounds or recent burns), thrombophlebitis and lack of sensitivity. Not contraindicated in the presence of implants

Patient sensations based on IAS value (IAS = INDIBA analog scale)

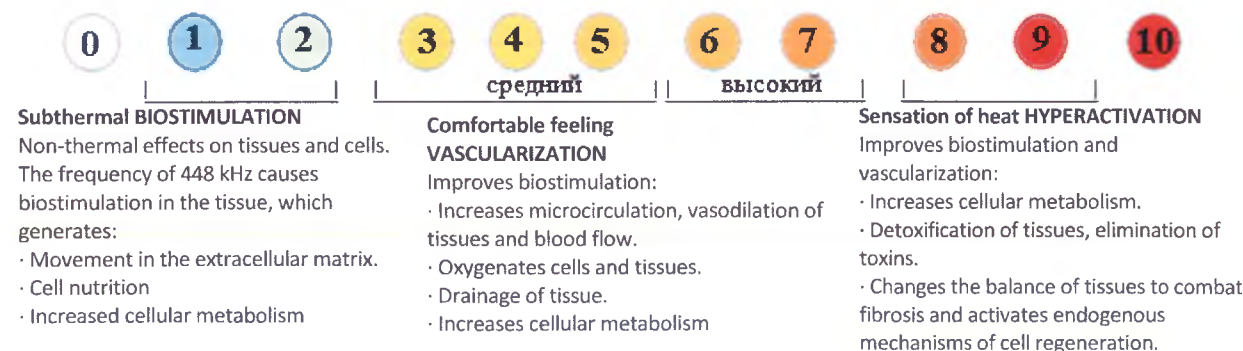
This protocol, which is exclusively owned by INDIBA, S.A., is protected by intellectual property law and may only be used by the rightful owners of the relevant equipment. Therefore, it is forbidden to reproduce, modify, distribute or transmit it in whole or in part without the written consent of INDIBA, S.A. Выпущенного в письменной форме. Все нарушения будут преследоваться по закону.

AFTER INTRODUCTION OF CALCIUM HYDROXYAPATITE, POLYCAPROLACTON OR POLYLACTIC ACID (FACE 20 min)

INDIBA
DEEP CARE
MEDICAL FREQUENT SYSTEM

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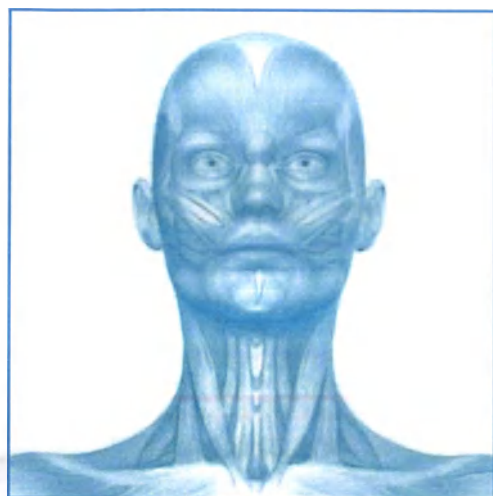POSITION

Lying on your back. Return flat electrode on the interscapular region (ensure good contact with the tissue)

DESCRIPTION OF STAGES OF THE PROCEDURE

CAP 10 min IAS* 6 – 8

RES 10 min IAS* 6 – 8

ELECTRODE CAP: 30 - 40 mm, RES 35 mm

APPLICATION DATA

Performing movements along the drainage line and towards the muscles

SESSIONS

Once a week

Total: 4 to 6 sessions

BEGINNING OF THE PROCEDURE

15-30 days after administration

NOTES

If necessary, a subthermal procedure (IAS 1 with signal modulation) can be performed immediately after the introduction of fillers to limit side effects from the introduction, without reducing the effectiveness of the procedure. This protocol will consist of 1 to 3 sessions with an interval of 24 hours between sessions in accordance with the needs of the patient.

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RECOMMENDATIONS

It is recommended to wait 30 days after the introduction of fillers to begin the procedure.

To achieve a long-term therapeutic effect, it is recommended to prescribe one course of treatment every 1 or 2 months.

Check for past hypersensitivity to face cream ingredients.

COMMENTS

The proposed protocols are general and should be adapted to the individual requirements or circumstances of the patient, taking into account age, hormone imbalance, concomitant conditions, lifestyle, quality of life and history of surgical intervention. All of these factors can affect the end result. For external use only.


APPLICATION TECHNOLOGY

- 1. Position:** hold the electrode evenly on the surface of the body; do not tilt.
- 2. Pressure:** maintain sufficient pressure to ensure adequate contact between the skin and the surface of the electrode. The pressure should be stable and not cause discomfort to the patient.
- 3. Speed:** move the electrode in a circular motion, each time slowly increasing the diameter of the circle. Movements should be slow (one circle per second) and smooth.
- 4. Температура:** ensure that the temperature is suitable for the patient and is comfortable.
- 5. Гигиена:** To optimize the transfer of energy to the patient when using the product, it is recommended that the user hold onto the electrode handle during the procedure. If it is necessary to use manual therapy methods, INDIBA representatives recommend the use of gloves (not included). This undoubtedly optimizes the transfer of energy to the patient and ensures the best result of the procedure.

CONTRAINDICATIONS

The procedure should not be carried out in the presence of a pacemaker or other electronic implants, in pregnant women, in case of skin wounds (open wounds or recent burns), thrombophlebitis and lack of sensitivity. Not contraindicated in the presence of implants

Patient sensations based on IAS value (IAS = INDIBA analog scale)

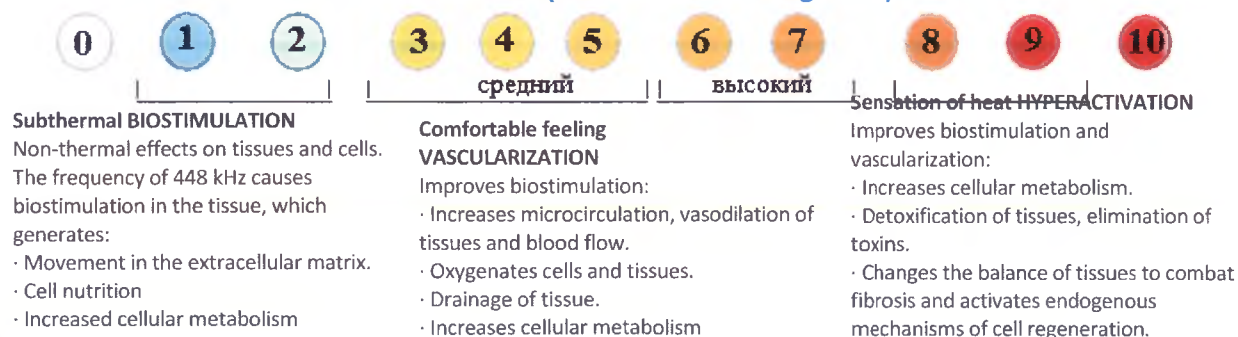
This protocol, which is exclusively owned by INDIBA, S.A., is protected by intellectual property law and may only be used by the rightful owners of the relevant equipment. Therefore, it is forbidden to reproduce, modify, distribute or transmit it in whole or in part without the written consent of INDIBA, S.A., issued in written form. All violations will be prosecuted by law.

BEFORE MESOTHERAPY WITH VITAMINS, PEPTIDES, NON-RETICULAR HYALURONIC ACID, PLASMOTHERAPY

FACE AND NECK(20 min)

INDIBA
DEEP CARE
MEDICAL PROIONIC SYSTEM

Ensuring the optimal state of tissue before the procedure, strengthening the appropriate intradermal distribution and assimilation of the infiltrated products due to the increased permeability of the cell membrane caused by the procedure. Therefore, the risk of hematomas due to infiltration will be reduced.

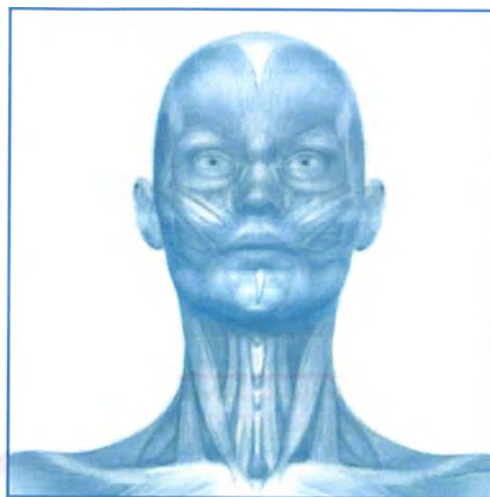POSITION

Lying on your back. Return flat electrode on the interscapular region (ensure good contact with the tissue)

DESCRIPTION OF STAGES OF THE PROCEDURE

CAP 10 min IAS* 6 – 8

RES 10 min IAS* 6 – 8

ELECTRODE CAP: 30 - 40 mm, RES 35 mm

APPLICATION DATA

Introduction 8 hours after the procedure to avoid the resulting thermal effect

SESSIONS

2-3 once a week

Total:2-4 sessions

BEGINNING OF THE PROCEDURE

1-2 weeks before administration

NOTES

If necessary, chemical peeling can be performed as an additional treatment 24/48 hours before the procedure.

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
AFTER MESOTHERAPY BY VITAMINS, PEPTIDES, NERETICULAR HYALURONIC ACID, PLASMOTHERAPY FACE AND NECK(20 min)

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Strengthening the biological effects of mesotherapy on the treated tissue (neocollagenesis, neoVASCULARIZATION, etc.) In addition, it is also possible to reduce the side effects of treatment (bruises, umoborizatsiya, edema ...).

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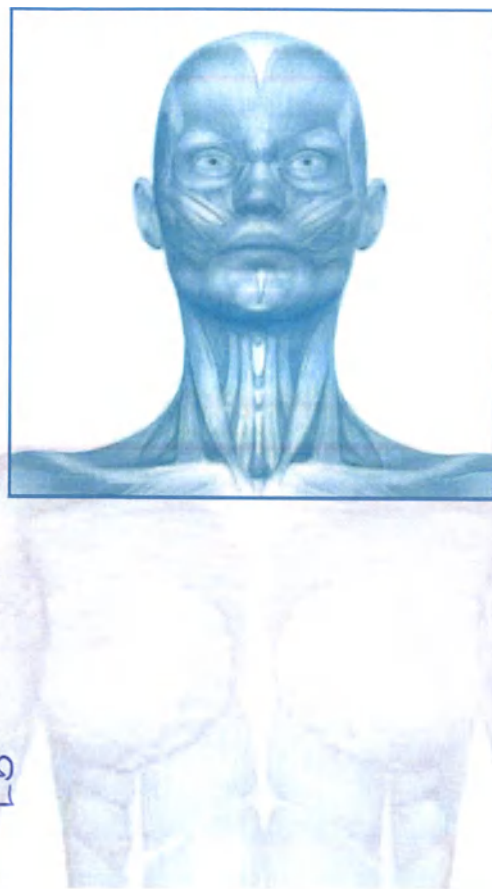15-30 days after administration

NOTES


In case of bruising, it is recommended to immediately stop the procedure, leaving at least 8 hours between sessions, depending on the needs of the patient. Perform circular motions in the area of bruises in CAP mode 5 - 7 min with IAS 1 - 2 with signal modulation.

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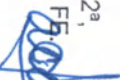AFTER THE INTRODUCTION OF BOTULINOTOXIN FOR REMOVING MIMIC WRINKLES FACE (20 min)

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An additional procedure to improve the condition of tissues. Promotes dermoepidermal regeneration after muscle relaxation, caused by the therapeutic effect of botulinum toxin with facial wrinkles. Improving the results.

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15 days after administration

NOTES

If necessary, chemical poring can be done 24-48 hours before the procedure.

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Sección 2ª

AFTER THE INTRODUCTION OF BOTULINOTOXIN FOR REMOVING MIMIC WRINKLES

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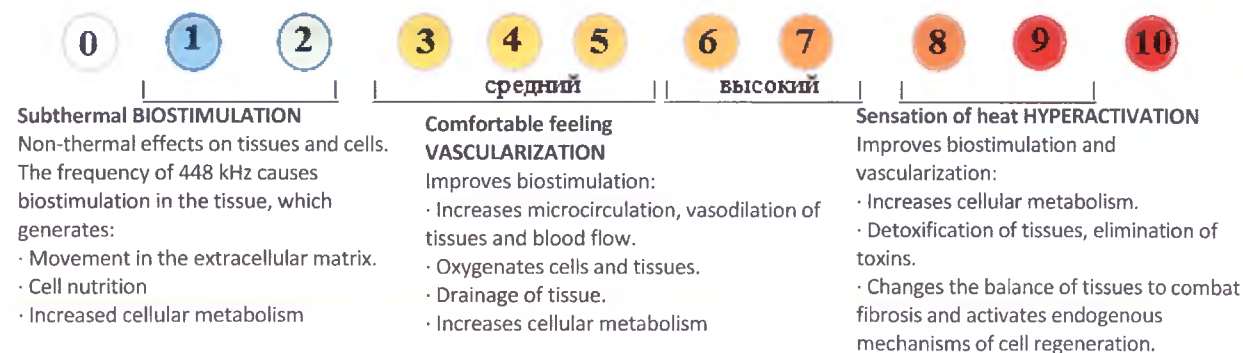
This protocol, which is exclusively owned by INDIBA, S.A., is protected by intellectual property law and may only be used by the rightful owners of the relevant equipment. Therefore, it is forbidden to reproduce, modify, distribute or transmit it in whole or in part without the written consent of INDIBA, S.A.. Any infringement will be prosecuted by law.

FOLLOWING TREATMENT OF SIDE EFFECTS CAUSED BY A NON-ABLATION FRACTION LASER AND FRACTIONAL ABLATION RADIO FREQUENCY

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FACE (20 min)

Control of the inflammatory and edematous process, without reducing its beneficial effects, accelerating the patient's recovery, improving results and reducing the risk of side effects (post-inflammatory hyperpigmentation ...).

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ELECTRODE RES 35 mm

APPLICATION DATA

Performing linear movements without shifting the tissue under a sterile hydrogel dressing and applying an electrically conductive cream under the electrode to slip

SESSIONS

1 once a week

Total: 4-6 sessions

BEGINNING OF THE PROCEDURE

15 days after administration

NOTES

If necessary, chemical poring can be done 24-48 hours before the procedure.

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
APPLICATION TECHNOLOGY

1. **Position:** hold the electrode evenly on the surface of the body; do not tilt.
2. **Pressure:** maintain sufficient pressure to ensure adequate contact between the skin and the surface of the electrode. The pressure should be stable and not cause discomfort to the patient.
3. **Speed:** move the electrode in a circular motion, each time slowly increasing the diameter of the circle. Movements should be slow (one circle per second) and smooth.
4. **Температура:** ensure that the temperature is suitable for the patient and is comfortable.
5. **Гигиена:** To optimize the transfer of energy to the patient when using the product, it is recommended that the user hold onto the electrode handle during the procedure. If it is necessary to use manual therapy methods, INDIBA representatives recommend the use of gloves (not included). This undoubtedly optimizes the transfer of energy to the patient and ensures the best result of the procedure. оцедурь.

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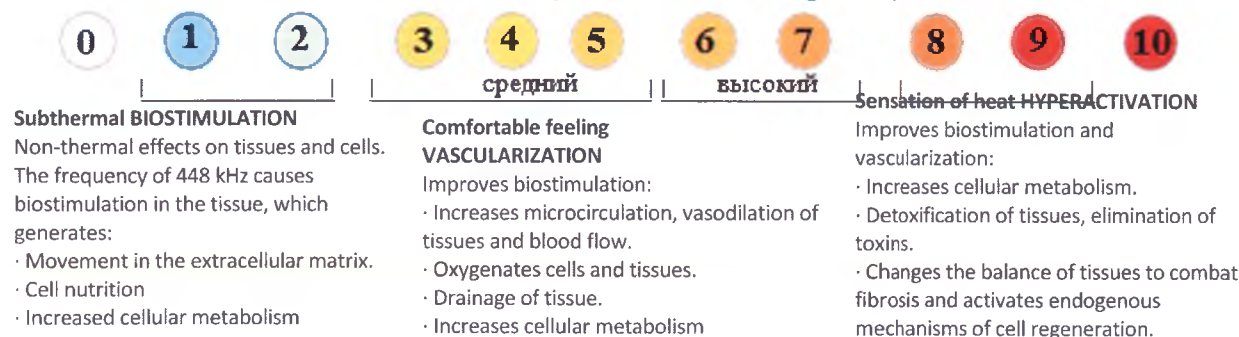This protocol, which is exclusively owned by INDIBA, S.A., is protected by intellectual property law and may only be used by the rightful owners of the relevant equipment. Therefore, it is forbidden to reproduce, modify, distribute or transmit it in whole or in part without the written consent of INDIBA, S.A.. Все нарушения будут преследоваться по закону.

NEXT TREATMENT OF SIDE EFFECTS CAUSED BY IFL AND A VASCULAR LASER (PDL, KTP, Nd: YAG, etc.)

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FACE (20 min)

Control of the inflammatory and edematous process, without reducing its beneficial effects, accelerating the patient's recovery, improving results and reducing the risk of side effects (post-inflammatory hyperpigmentation ...).

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ELECTRODE CAP 30-40 mm, RES 35 mm

APPLICATION DATA

Perform linear or circular motions without slipping over fabric

SESSIONS

2-3 once a week

Total: 1-3 sessions

BEGINNING OF THE PROCEDURE

24 hours after exposure

NOTES

Immediate treatment after laser exposure is recommended, leaving at least 8 hours between sessions according to patient needs

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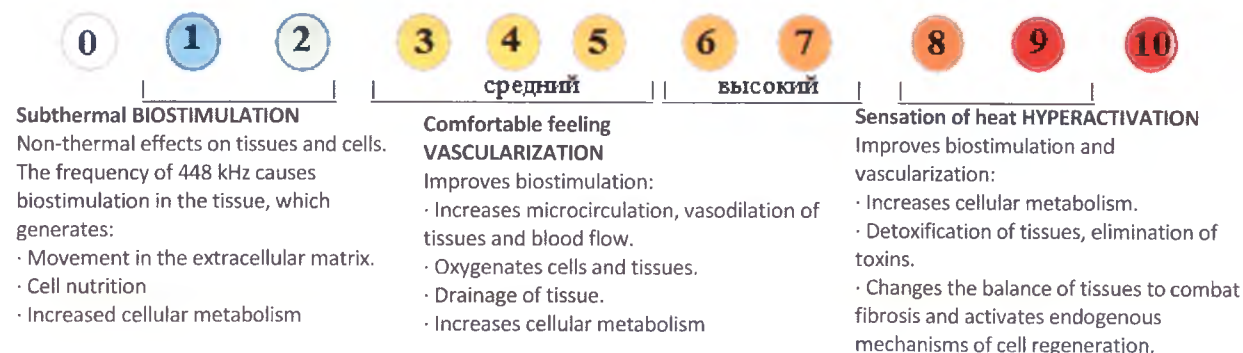
This protocol, which is exclusively owned by INDIBA, S.A., is protected by intellectual property law and may only be used by the rightful owners of the relevant equipment. Therefore, it is forbidden to reproduce, modify, distribute or transmit it in whole or in part, to copy, to disseminate or to use it for any other purpose without the prior written consent of INDIBA, S.A. All violations will be pursued by law.

BEFORE MESOTHERAPY AND LASER PROCEDURE (DEPIGMENTATION, NON-ABLATION ETC.) IN THE INFRAORBITAL AREA (SACKS UNDER EYES)

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INFRAORBITAL AREA (10 min)

Facilitate drainage of the area to be infiltrated to avoid bruising and swelling of the eyelids.

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24-48 hours after exposure

NOTES

If necessary, a subthermal follow-up procedure (IAS 1 with signal modulation) can be performed immediately after treatment if a hematoma occurs.

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
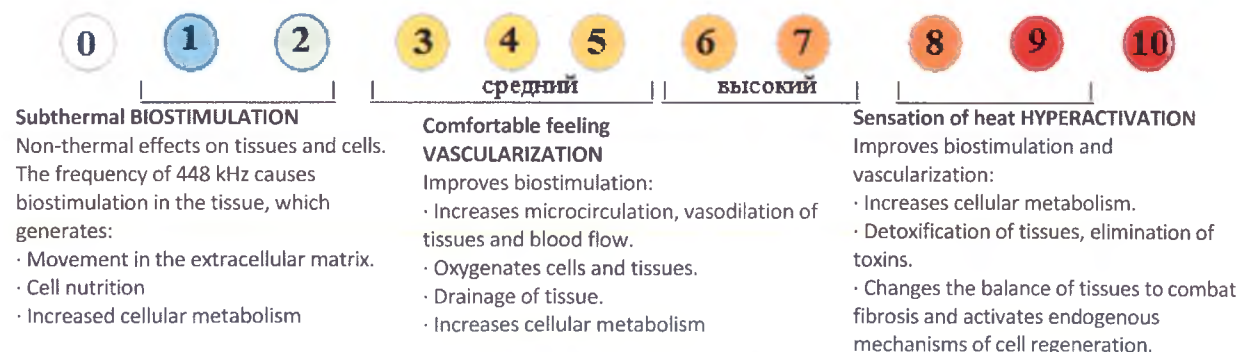
This protocol, which is exclusively owned by INDIBA, S.A., is protected by intellectual property law and may only be used by the rightful owners of the relevant equipment. Therefore, it is forbidden to reproduce, modify, distribute or transmit it in whole or in part without the written consent of INDIBA, S.A.. Все нарушения будут преследоваться по закону.

BEFORE CAPILLARY SEALING USING VARIOUS TECHNOLOGIES (MESOTHERAPY, CARBOXYTHERAPY, TOP, PHOTO REJUVENATION ...)

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Fronto-dark and occipital region (15 min.)

Immediate tissue compaction, eliminating product non-dispersion.

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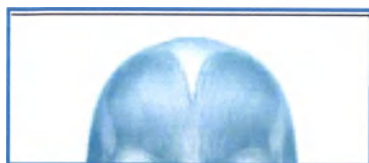Consider the possible decrease in the effectiveness of the procedure in patients with collagen diseases or in patients taking high doses of corticosteroids.

Check for past hypersensitivity to face cream ingredients Face Cream.

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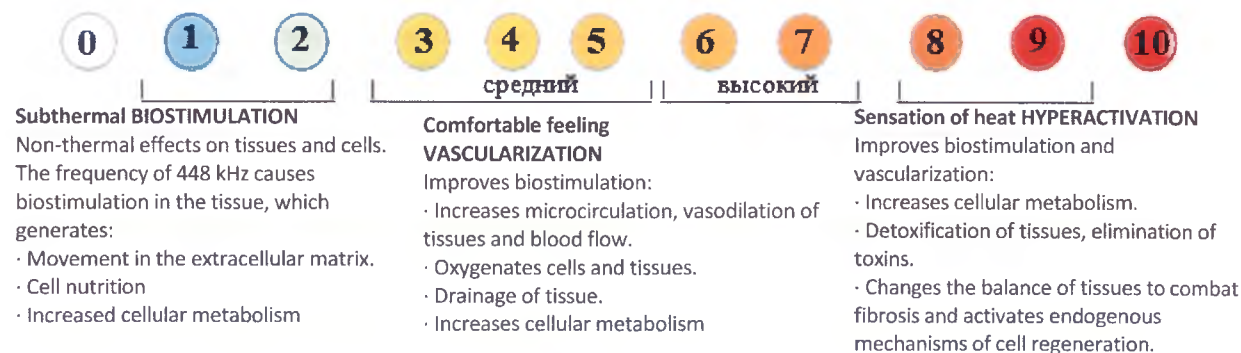
This protocol, which is exclusively owned by INDIBA, S.A., is protected by intellectual property law and may only be used by the rightful owners of the relevant equipment. Therefore, it is forbidden to reproduce, modify, distribute or transmit it in whole or in part without the prior written consent of INDIBA, S.A.. Выпущено в письменной форме. Все нарушения будут преследоваться по закону.

AFTER CAPILLARY SEALING USING VARIOUS TECHNOLOGIES (MESOTHERAPY, CARBOXYTHERAPY, TOP, PHOTO REJUVENATION ...)

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Fronto-dark and occipital region (15 min.)

Control of inflammatory and edematous processes, without reducing beneficial effects.

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72 hours after exposure

NOTES

It is recommended to start the protocol 72 hours after treatment, leaving at least 8 hours between sessions according to the needs of the patient

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RECOMMENDATIONS

After 1 month after treatment, INDIBA procedures can be used as an adjuvant (IAS 3-5) against aging.

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
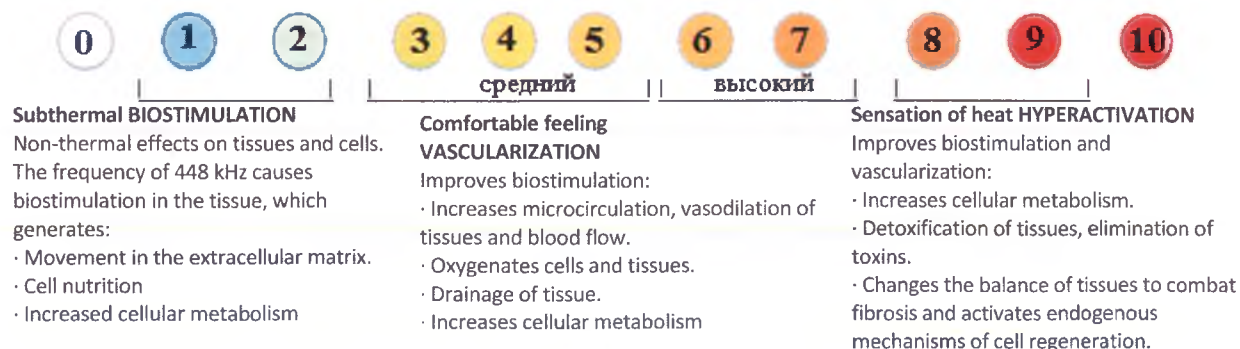
This protocol, which is exclusively owned by INDIBA, S.A., is protected by intellectual property law and may only be used by the rightful owners of the relevant equipment. Therefore, it is forbidden to reproduce, modify, distribute or transmit it in whole or in part without the written consent of INDIBA, S.A., expressed in writing. All violations will be pursued by law.

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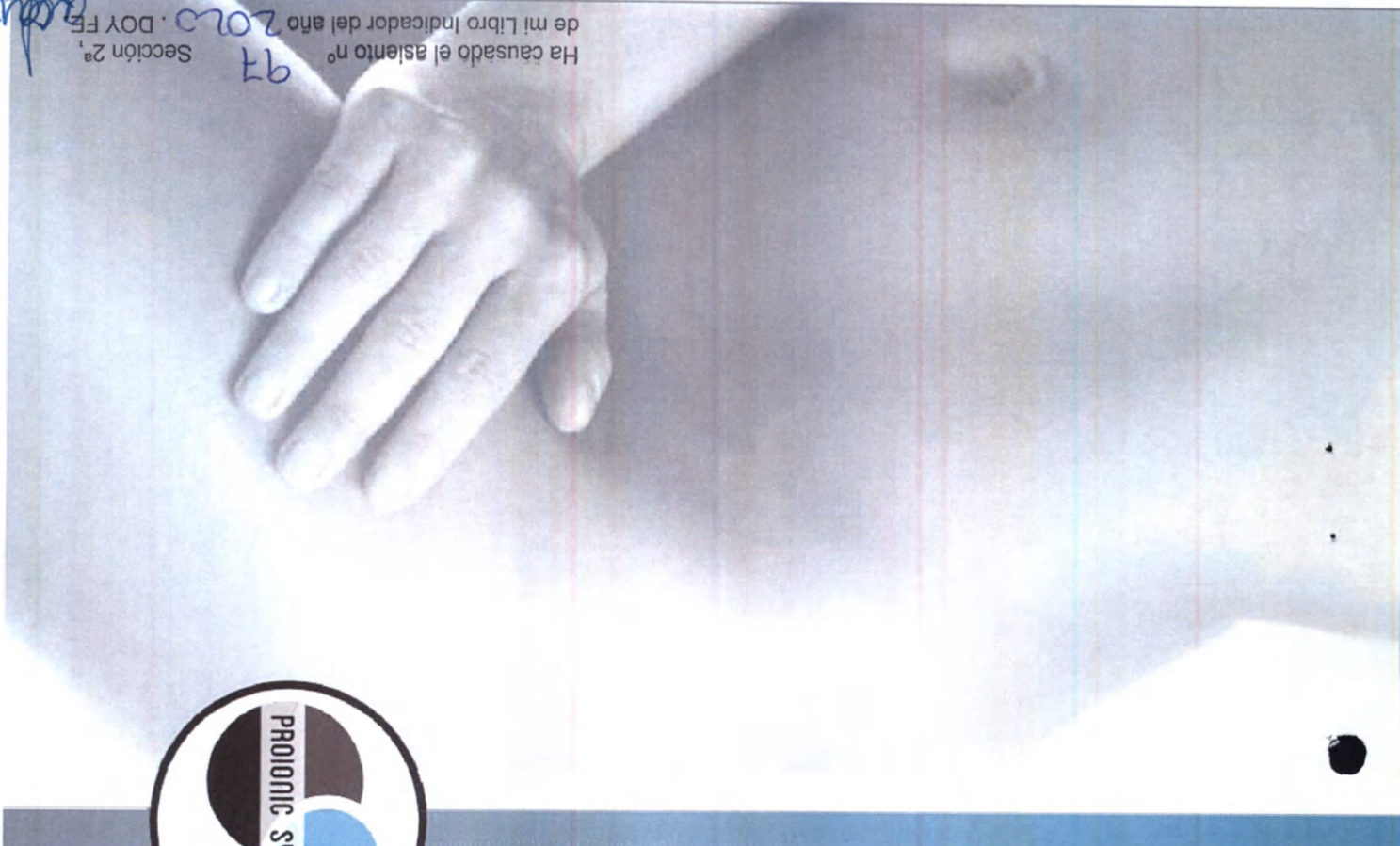TRN0239 V2 2018.06.26 EN

GENITAL REJUVENATION (20 min.)

Strengthening, enhancing hydration, closing the vagina, tightening collagen

POSITION

Lithotomy position. Reversible flat electrode on the lumbosacral

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RES 35 mm: 5 min IAS* 8 from each side

APPLICATION DATA

Perform circular or linear movements depending on the patient's anatomy.

First process one half of the vaginal opening and then the other. The temperature should be around 40 - 42 ° C.

SESSIONS

1 time per month

Total: 3 sessions

NOTES

ELECTRODE RES, ELECTRODE RES Longitudinal: for individual use only. Clean ELECTRODE twice with soapy water. Then apply a disinfectant that is intended for professional use (phenol or glutaraldehyde 2%) for 30 min, and then rinse with water.

ELECTRODE CAP, ELECTRODE CAP longitudinal: for individual use only. Clean ELECTRODE twice with soapy water.

Disinfect the handles for ELECTRODES and returnable ELECTRODES using the products used in the hospital. If in doubt, contact the manufacturer.

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RECOMMENDATIONS

We recommend this procedure after more targeted therapy or as a concomitant treatment to improve the end result. Treatment should be accompanied by exercise and a healthy diet.

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
APPLICATION TECHNOLOGY

1. **Position:** hold the electrode evenly on the surface of the body; do not tilt.
2. **Pressure:** maintain sufficient pressure to ensure adequate contact between the skin and the surface of the electrode. The pressure should be stable and not cause discomfort to the patient.
3. **Speed:** move the electrode in a circular motion, each time slowly increasing the diameter of the circle. Movements should be slow (one circle per second) and smooth.
4. **Temperature:** ensure that the temperature is suitable for the patient and is comfortable.
5. **Hygiene:** INDIBA representatives recommend the use of nitrile or latex disposable gloves (not supplied).

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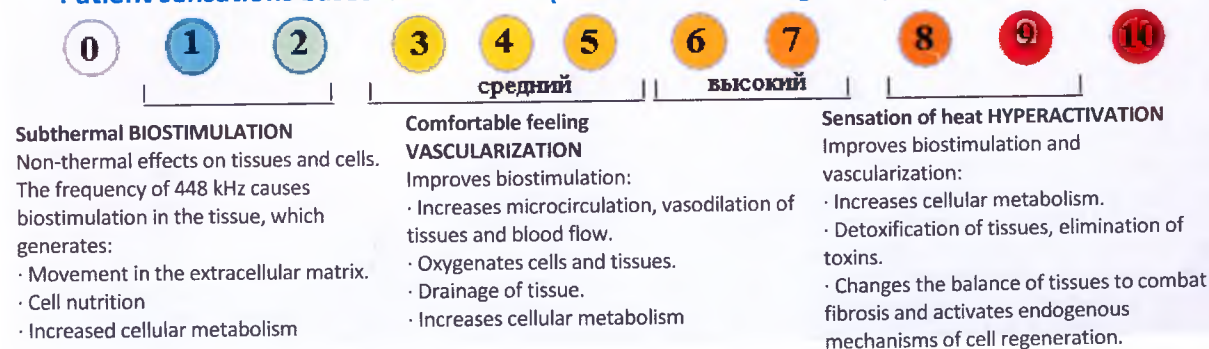which is exclusively owned by INDIBA, S.A., is protected by intellectual property law and may only be used by the rightful owners of the relevant equipment. Therefore, it is forbidden to reproduce, modify, distribute or transmit it in whole or in part without the written consent of INDIBA, S.A.. Выпускено в письменной форме. Все нарушения будут преследоваться по закону.

BEFORE GENITAL SURGERY (20 min)

Preparation of tissue before surgery, increased blood supply and oxygenation.

POSITION

Lithotomy position. Reversible flat electrode on the lumbosacral

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RES 35 mm: 5 min IAS* 6-7 from each side

APPLICATION DATA

Perform circular or linear movements depending on the patient's anatomy.

First process one half of the vaginal opening and then the other.

SESSIONS

Total: 1 or 2 sessions

The last session should be performed 24-48 hours before surgery.

NOTES

ELECTRODE RES, ELECTRODE RES Longitudinal: for individual use only. Clean ELECTRODE twice with soapy water. Then apply a disinfectant that is intended for professional use (phenol or glutaraldehyde 2%) for 30 min, and then rinse with water.

ELECTRODE CAP, ELECTRODE CAP longitudinal: for individual use only. Clean ELECTRODE twice with soapy water.

Disinfect the handles for ELECTRODES and returnable ELECTRODES using the products used in the hospital. If in doubt, contact the manufacturer.

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RECOMMENDATIONS

We recommend this procedure after more targeted therapy or as a concomitant treatment to improve the end result. Treatment should be accompanied by exercise and a healthy diet.

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
APPLICATION TECHNOLOGY

1. **Position:** hold the electrode evenly on the surface of the body; do not tilt.
2. **Pressure:** maintain sufficient pressure to ensure adequate contact between the skin and the surface of the electrode. The pressure should be stable and not cause discomfort to the patient.
3. **Speed:** move the electrode in a circular motion, each time slowly increasing the diameter of the circle. Movements should be slow (one circle per second) and smooth.
4. **Temperature:** ensure that the temperature is suitable for the patient and is comfortable.
5. **Hygiene:** INDIBA representatives recommend the use of nitrile or latex disposable gloves (not supplied).

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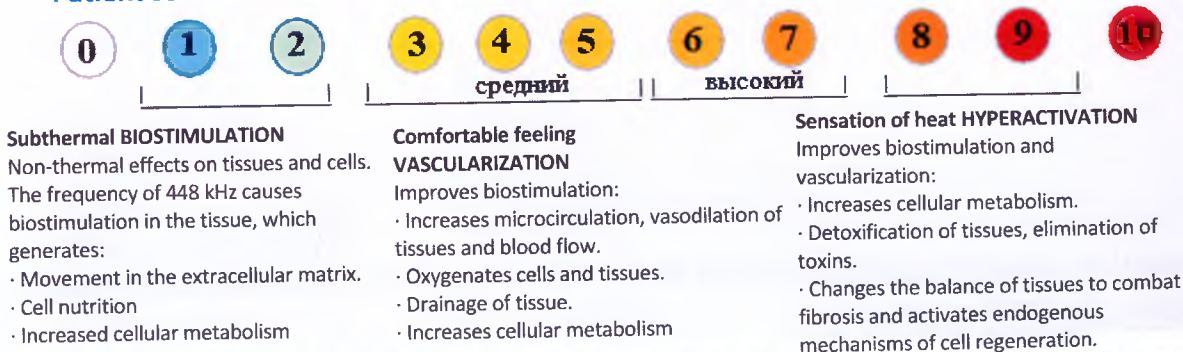which is exclusively owned by INDIBA, S.A., is protected by intellectual property law and may only be used by the rightful owners of the relevant equipment. Therefore, it is forbidden to reproduce, modify, distribute or transmit it in whole or in part without the written consent of INDIBA, S.A. Все нарушения будут преследоваться по закону.

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Accelerate wound healing, prevent adhesions and fibrosis, strengthen tissues and improve hydration.

POSITION

Lithotomy position. Reversible flat electrode on the lumbosacral

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CAP 40 mm: 2 min IAS* 1 from each side

APPLICATION DATA

Perform circular or linear movements depending on the patient's anatomy. Perform the procedure only after the inflammatory process disappears. This protocol can be continued with the genital rejuvenation protocol. If necessary, you can combine this protocol with others.

SESSIONS

Once a week

Total: 4 sessions

NOTES

ELECTRODE RES, ELECTRODE RES Longitudinal: for individual use only.

Clean ELECTRODE twice with soapy water. Then apply a disinfectant that is intended for professional use (phenol or glutaraldehyde 2%) for 30 min, and then rinse with water.

ELECTRODE CAP, ELECTRODE CAP longitudinal: for individual use only.

Clean ELECTRODE twice with soapy water.

Disinfect the handles for ELECTRODES and returnable ELECTRODES using the products used in the hospital. If in doubt, contact the manufacturer.

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RECOMMENDATIONS

We recommend this procedure after more targeted therapy or as a concomitant treatment to improve the end result. Treatment should be accompanied by exercise and a healthy diet.

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
APPLICATION TECHNOLOGY

1. **Position:** hold the electrode evenly on the surface of the body; do not tilt.
2. **Pressure:** maintain sufficient pressure to ensure adequate contact between the skin and the surface of the electrode. The pressure should be stable and not cause discomfort to the patient.
3. **Speed:** move the electrode in a circular motion, each time slowly increasing the diameter of the circle. Movements should be slow (one circle per second) and smooth.
4. **Temperature:** ensure that the temperature is suitable for the patient and is comfortable.
5. **Hygiene:** INDIBA representatives recommend the use of nitrile or latex disposable gloves (not included).

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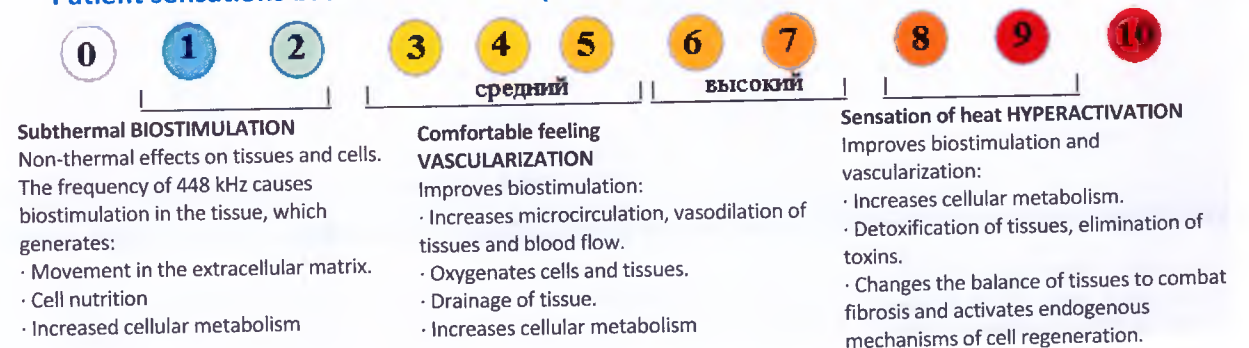
which is exclusively owned by INDIBA, S.A., is protected by intellectual property law and may only be used by the rightful owners of the relevant equipment. Therefore, it is forbidden to reproduce, modify, distribute or transmit it in whole or in part, or to use it for any other purpose without the prior written consent of INDIBA, S.A. Выпускенного в письменной форме. Все нарушения будут преследоваться по закону.

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Accelerate wound healing, prevent adhesions and fibrosis, strengthen tissues and improve hydration.

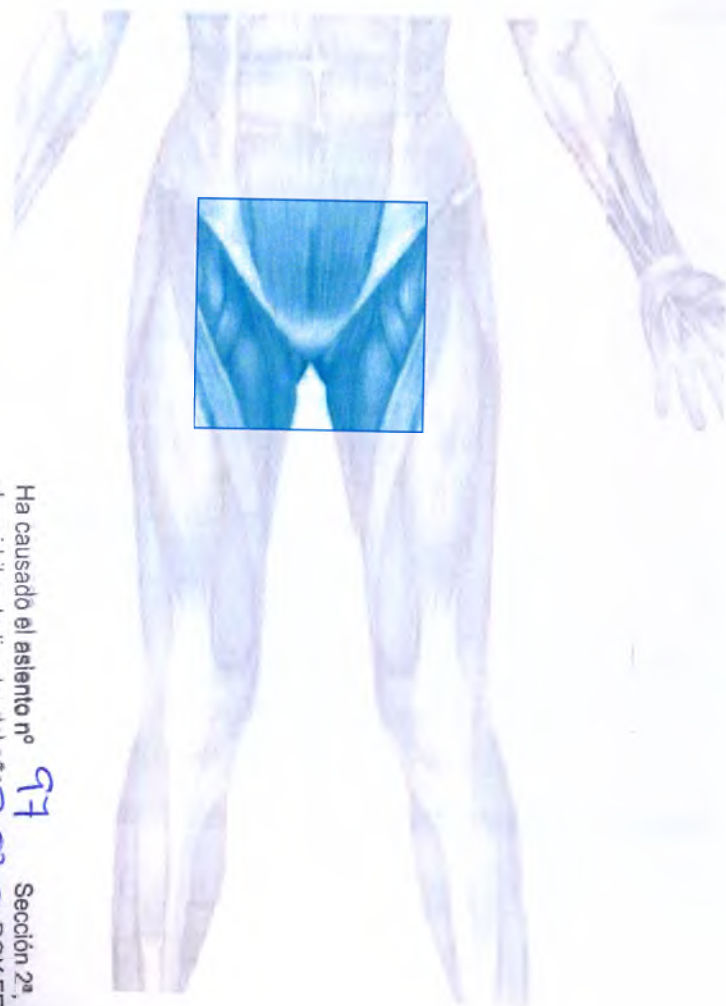POSITION

Lithotomy position. Reversible flat electrode on the lumbosacral

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CAP 40 mm: 2 min IAS* 1 from each side

APPLICATION DATA

Perform circular or linear movements depending on the patient's anatomy. Perform the procedure only after the inflammatory process disappears. This protocol can be continued with the genital rejuvenation protocol. If necessary, you can combine this protocol with others.

SESSIONS

Once a week

Total: 4 sessions

NOTES

ELECTRODE RES, ELECTRODE RES Longitudinal: for individual use only.

Clean ELECTRODE twice with soapy water. Then apply a disinfectant that is intended for professional use (phenol or glutaraldehyde 2%) for 30 min, and then rinse with water.

ELECTRODE CAP, ELECTRODE CAP longitudinal: for individual use only.

Clean ELECTRODE twice with soapy water.

Disinfect the handles for ELECTRODES and returnable ELECTRODES using the products used in the hospital. If in doubt, contact the manufacturer.

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INDIBA
DEEP CARI
MEDICAL FREONIC SYSTEM

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We recommend this procedure after more targeted therapy or as a concomitant treatment to improve the end result. Treatment should be accompanied by exercise and a healthy diet.

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The proposed protocols are general and should be adapted to the individual requirements or circumstances of the patient, taking into account age, hormone imbalance, concomitant conditions, lifestyle, quality of life and history of surgical intervention. All of these factors can affect the end result.

For external use only.

1. Position: hold the electrode evenly on the surface of the body; do not tilt.
2. Pressure: maintain sufficient pressure to ensure adequate contact between the skin and the surface of the electrode. The pressure should be stable and not cause discomfort to the patient.
3. Speed: move the electrode in a circular motion, each time slowly increasing the diameter of the circle. Movements should be slow (one circle per second) and smooth.
4. Temperature: ensure that the temperature is suitable for the patient and is comfortable.
5. Hygiene: INDIBA representatives recommend the use of nitrile or latex disposable gloves (not supplied).

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Non-thermal effects on tissues and cells.
The frequency of 448 kHz causes biostimulation in the tissue, which generates:

- Movement in the extracellular matrix.
- Cell nutrition
- Increased cellular metabolism

- Increases microcirculation, vasodilation of tissues and blood flow.
- Oxygenates cells and tissues.
- Drainage of tissue.
- Increases cellular metabolism

- Improves biostimulation and vascularization:
- Increases cellular metabolism.
- Detoxification of tissues, elimination of toxins.
- Changes the balance of tissues to combat fibrosis and activates endogenous mechanisms of cell regeneration.

This protocol, which is exclusively owned by INDIBA, S.A., is protected by intellectual property law and may only be used by the rightful owners of the relevant equipment. Therefore, it is forbidden to reproduce, modify, distribute or transmit it in whole or in part without the prior written consent of INDIBA, S.A.

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
Reduced sagging tissues, improved hydration and the elimination of dryness, rejuvenation of epithelial tissue and vaginal mucosa. Decreased vaginal atrophy due to premenstrual syndrome. Improving local estrogen trophism. Increased elasticity and firmness of the vagina. Activation of neocollagenesis. Stimulation of contractility and thickening of the walls of the vagina.

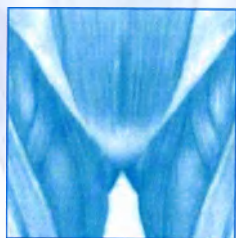POSITION

Lithotomy position. Reversible flat electrode on the lumbosacral

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External application: CAP 30-40 mm: 2 min IAS*8 from each side

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APPLICATION DATA

External use: Perform circular or linear movements depending on the patient's anatomy. Maintain sufficient pressure to ensure perfect contact between the skin and the entire surface of ELECTRODE.

Internal use: Use a 1-minute interval with pressure on the sides of the vagina (upper, right, left, lower): a total of 3 cycles of 4 minutes (1 minute from each side) to complete the procedure. ELECTRODE should be fully in contact with the vaginal wall. Before starting the procedure, ELECTRODE must be completely immersed in the vagina. Be sure to stop the procedure before removing the ELECTRODE.

SESSIONS

Once a week

Total: 4 sessions

NOTES

ELECTRODE RES, ELECTRODE RES Longitudinal: for individual use only. Clean ELECTRODE twice with soapy water. Then apply a disinfectant that is intended for professional use (phenol or glutaraldehyde 2%) for 30 min, and then rinse with water.

ELECTRODE CAP, ELECTRODE CAP longitudinal: for individual use only. Clean ELECTRODE twice with soapy water.

Disinfect the handles for ELECTRODES and returnable ELECTRODES using the products used in the hospital. If in doubt, contact the manufacturer.

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97

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[Signature]

AUXILIARY TREATMENT OF URINE STRESS CONTINUITY

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
RECOMMENDATIONS

We recommend this procedure after more targeted therapy or as a concomitant treatment to improve the end result. Treatment should be accompanied by exercise and a healthy diet.

APPLICATION TECHNOLOGY

1. **Position:** hold the electrode evenly on the surface of the body; do not tilt.
2. **Pressure:** maintain sufficient pressure to ensure adequate contact between the skin and the surface of the electrode. The pressure should be stable and not cause discomfort to the patient.
3. **Speed:** move the electrode in a circular motion, each time slowly increasing the diameter of the circle. Movements should be slow (one circle per second) and smooth.
4. **Temperature:** ensure that the temperature is suitable for the patient and is comfortable.
5. **Hygiene:** INDIBA representatives recommend the use of nitrile or latex disposable gloves (not supplied).

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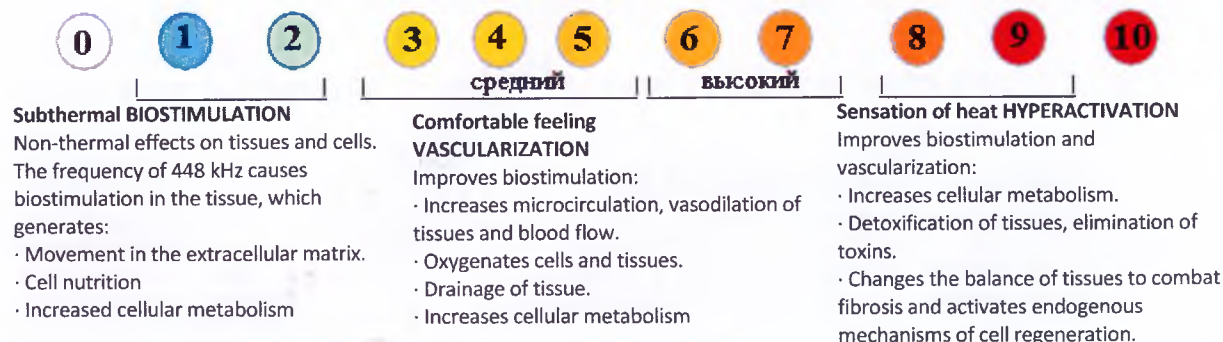
PROTOCOL FOR THE TREATMENT OF POST-PERMANENT INJURY (NOT IN ACUTE PHASE) (20min)

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Speed up the healing process. Fibrosis treatment. Pain relief. Restoring vaginal elasticity. Activation of neocalogenesis. Stimulation of contractility and activation of local neocollagenesis. Improving tissue elasticity.

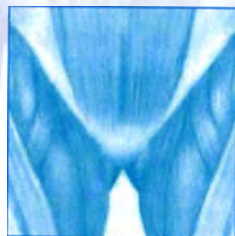POSITION

Lithotomy position. Reversible flat electrode on the lumbosacral

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External application: CAP 30-40 mm: 5 min IAS* 6-7

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APPLICATION DATA

External use: Perform circular or linear movements depending on the patient's anatomy. Maintain sufficient pressure to ensure perfect contact between the skin and the entire surface of ELECTRODE.

Internal application: ELECTRODE should be fully in contact with the vaginal wall. Before starting the procedure, ELECTRODE must be completely immersed in the vagina. Be sure to stop the procedure before removing the ELECTRODE.

SESSIONS

Once a week

Total: 8 sessions

Other specific treatments are also recommended.

NOTES

ELECTRODE RES, ELECTRODE RES Longitudinal: for individual use only. Clean ELECTRODE twice with soapy water. Then apply a disinfectant that is intended for professional use (phenol or glutaraldehyde 2%) for 30 min, and then rinse with water.

ELECTRODE CAP, ELECTRODE CAP longitudinal: for individual use only. Clean ELECTRODE twice with soapy water.

Disinfect the handles for ELECTRODES and returnable ELECTRODES using the products used in the hospital. If in doubt, contact the manufacturer.

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PROTOCOL FOR THE TREATMENT OF POST-PERMANENT INJURY (NOT IN ACUTE PHASE)

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Alfonso Carbonell

This protocol, which is exclusively owned by INDIBA, S.A., is protected by intellectual property law and may only be used by the rightful owners of the relevant equipment. Therefore, it is forbidden to reproduce, modify, distribute or publish it in any form or by any means, electronic or mechanical, including photocopying, recording, or by any information storage and retrieval system. All rights reserved. INDIBA, S.A., выпускаемого в письменной форме. Все нарушения будут преследоваться по закону.

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This Apostille only certifies the authenticity of the signature and the capacity of the person who has signed the public document, and, where appropriate, the identity of the seal or stamp which the public document bears.

This Apostille does not certify the content of the document for which it was issued.

[This Apostille is not valid for use anywhere within Spain]

[To verify the issuance of this Apostille, see <https://sede.mjusticia.gob.es/eregister>]

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Copyright 2010 Indiba E. de. All rights reserved. Made in Spain.
Copyright© 2015 Indiba S.A. All rights reserved. Made in Spain