PROVINCE OF ONTARIO
Regional Municipality
of Waterloo

TO WIT

29623
TO ALL WHOM THESE PRESENTS
MAY COME, BE SEEN OR KNOWN


I, Victor Séan Gomes, a Notary Public, in and for the Province of Ontario, by Royal Authority duly appointed, residing at the City of Guelph, in the said Province,

DO CERTIFY AND ATTEST that the paper-writing hereto annexed is an original document purporting to be

a Unitron Letter, dated December 11, 2017 and signed by Amir Elahi-Panah, Manager, Quality Assurance, with attached Appendices 1 to 6, inclusive

an act whereof being requested I have granted under my Notarial Form and Seal of Office to serve and avail as occasion shall or may require.

IN TESTIMONY WHEREOF I have hereto subscribed my name and affixed my Notarial Seal of Office at Kitchener, Ontario, in the Province of Ontario, this 13th day of December, 2017.

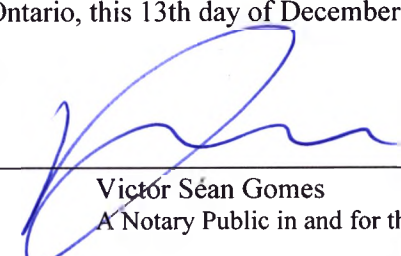Victor Séan Gomes
A Notary Public in and for the Province of Ontario

Date: December 11, 2017

To whom it may concern:

The following documents are provided for homologations purposes:

1. Appendix 1 – 029-5905-02f_uMic_Gd_ENG
2. Appendix 2 – 029-5989-02b_RemoteControl2_Gd_ENG
3. Appendix 3 – 029-5990-02b_uDirect3_Gd_ENG
4. Appendix 4 – 029-5992-02d_uTV3_Gd_ENG
5. Appendix 5 – 029-5993-02c_uStream_Gd_ENG
6. Appendix 6 – english-completeguide-icubeii

If you require further explanations or clarifications please contact the undersigned at Amir.Panah@Unitron.com


Best regards,

Amir Elahi-Panah
Manager, Quality Assurance

Appendix 1 – 029-5905-02f_uMic_Gd_ENG

uMic
guide

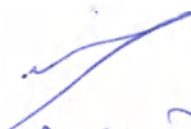Amir Panah

unitron. Hearing matters

Thank you

Thank you for choosing the Unitron uMic™

At Unitron, we care deeply about people with hearing loss. We work closely with hearing healthcare professionals to make advanced, purpose-driven solutions available to everyone. Because hearing matters.

Your uMic

Hearing healthcare professional: _____

Telephone: _____

Serial number of uMic: _____

Warranty: _____


Date of purchase: _____

Quick reference

Charging uMic


On/Off

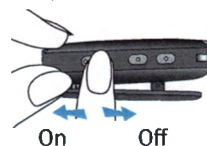Bluetooth pairing

Wearing uMic

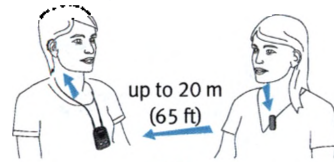Indicator


Switching on		Green 2 seconds
Switching off		Red 2 seconds
Charging		Solid red
Low battery warning		Additional red short blink as well as the main light
Streaming		Solid blue
Streaming paused		Short blink blue
Pairing mode		Rapid blink blue

Table of contents

Your uMic at a glance	2
Getting started	4
Charging the uMic	5
Turning the uMic on and off.....	7
Bluetooth pairing with the uMic.....	8
Understanding the status light	10
Functional test of uMic	12
Wearing the uMic correctly.....	13
Setting the uMic volume	14
Switching the uMic off.....	15
uMic and phone calls	16
uMic and watching TV	17
uMic streaming distance.....	18
Deleting all Bluetooth pairings.....	19
Use of the lanyard	20
Troubleshooting guide	21
Warnings	26
Information and explanation of symbols.....	31
Compliance information.....	33

Your uMic at a glance

- 1** Volume up (+)
- 2** Volume down (−)
- 3** Power switch (U/I)
- 4** Charging input (mini-USB)
- 5** Microphone opening
- 6** Status indicator light
- 7** Eyelet for lanyard
- 8 Retaining clip**

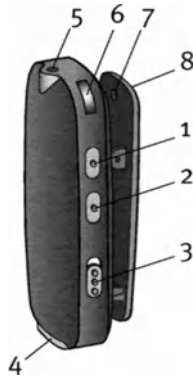Accessories

(Contents may vary depending on country)

Dual mini-USB splitter (1 in / 2 out)

Lanyard

Protective pouch

Other optional accessories

Unitron mini-USB power supply

12V car charger adapter

 uMic requires a uDirect 2, uDirect 3 or uStream, in order to be used with your Unitron hearing aid. See the appropriate user guide for details.

Getting started


The system consists of three main elements:

- uMic
- Unitron hearing aids
- A streamer

The uMic is compatible with three streaming devices, uDirect 2, uDirect 3 and uStream.

The uMic transmits the speaker's voice wirelessly to your hearing aids via the streamer.

The hearing aid microphone level will automatically be adjusted while the uMic signal is received. The level of adjustment is set by your hearing healthcare professional.

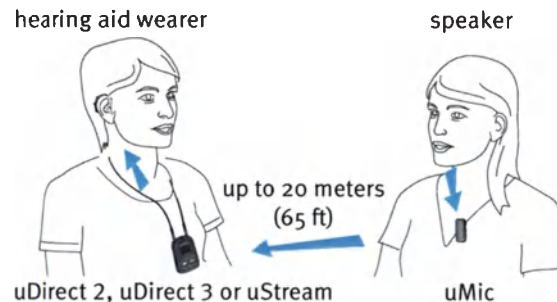4

Charging the uMic

Your uMic has a built-in, non-removable rechargeable battery. When the battery is low, the indicator light will blink red.

To charge:

1. Insert the Unitron power supply provided with the uDirect 2 or uDirect 3 into the uMic mini-USB socket.

2. Plug the power supply into a power outlet.

- Typical charging time is two hours. Charging is complete when the indicator light turns solid green.
- When charging your uMic for the first time, let it charge for at least three hours, even if the battery indicator is not lit or turns green earlier.

5

Using one power supply for two devices:

Insert the dual mini-USB splitter between the power supply and the devices as shown below.

With uDirect 2, uDirect 3 or uStream:

Both devices can be charged at the same time.

With uTV™ 3:


The uMic can be charged while uTV 3 is being used.

 It is safe to leave the uMic connected to the charger overnight.

Turning the uMic on and off

Slide the power switch  as shown to turn the uMic on or off.

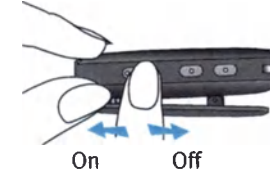Sleep mode:

If after 30 minutes no streamer device is found the uMic will power down to save battery life even if the power switch is in the on position.

To recover from sleep mode, switch the uMic off and turn it on again.

Bluetooth pairing with the uMic

Put uMic into pairing mode:

 If the uMic has never been paired before, pairing mode will be automatically activated after turning it on.

To manually start pairing mode:

1. Switch the uMic off.
2. Press and hold the volume up button (+).
3. Turn the uMic on while still pushing the (+) button.
4. Release the (+) button after two seconds when the indicator light rapidly blinks blue.
5. Make sure the streamer and uMic are within 1 m (3 ft) range of each other.

The uMic will remain in pairing mode for two minutes or until the pairing is completed.

Completion of pairing:

Wait until both indicator lights turn solid blue. This can take up to two minutes but typically is done within a few seconds.

-  The uMic can be paired to a maximum of four devices. A fifth pairing will replace the first paired device.
-  Two devices need to be paired only once. They remain paired after switching off and will automatically connect after being turned on again.

Understanding the status light

Information during operation:

Switching on	Green 2 seconds	
Streaming	Solid blue	
Streaming paused	Short blue blink	
Pairing mode	Rapid blue blink	

Information about power and battery:

Switching on	Green 2 seconds	
Switching off	Red 2 seconds	
Charging	Solid red	
Charged and full battery	Solid green	
Low battery warning	Additional red short blink as well as the main light	
uMic sleep mode or off	Indicator is unlit	

Switch the uMic off and on again to turn the device on from sleep mode.

If the red “low battery” blinks:

- Battery level is <20%
- Less than 1.5 hours are left for streaming
- Consider charging the uMic

The uMic can also be charged during operation.

Functional test of uMic

After pairing to your streamer, perform a quick functional test by gently tapping on the microphone opening.


You should clearly hear the tapping noise in your hearing aids.

 Line of sight between the uMic and the streamer is usually not required. If the audio transmission ends abruptly, refer to the troubleshooting guide.

Wearing the uMic correctly

1. Gently press to open the clip.
2. Ensure the uMic is clipped securely to clothing as shown in the picture.
3. To prevent scratching sounds, do not cover the uMic with clothing during use.

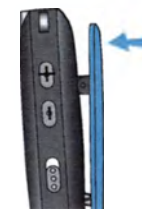Distance to the speaker's mouth:

- Best results are achieved at 20 cm (8")
- Do not exceed 50 cm (20") for consistent audio

 The higher the ambient noise, the closer the uMic should be positioned to the speaker's mouth.

 Do not use a neckloop lanyard during streaming. Scratching noise may negatively impact the sound.

Setting the uMic volume

To change the uMic volume:

1. Speak into the uMic from a 20 cm (8") distance.
2. Use the volume (+/−) buttons to set the uMic volume to adjust the sensitivity and strength of the microphones.

If your hearing aids pick up too much ambient noise, reduce the volume on your hearing aids and increase the volume on your uMic.

-  The uMic volume setting will be stored even when switched off.
-  Your hearing healthcare professional can change the level of the hearing aid microphones while using the uMic to adjust the volume of sounds around you.

Switching the uMic off

To end the uMic transmission completely, switch off the uMic by sliding the power switch to the off position.

-  In the off position, the white marking cannot be seen.

Switch the uMic off when not in use.

-  If the uMic connection to a streamer was lost for more than 30 minutes, the uMic will go into sleep mode to save power. The light indicator will turn off and a uMic re-start (turn off and on again) is required to resume the transmission.

uMic and phone calls

Phone calls via Bluetooth will always take priority over the uMic streaming.

If uDirect 2, uDirect 3 or uStream streamer is paired with your mobile phone, you can still receive phone calls, even when using uMic. In case of an incoming call, the streamer automatically stops transmitting the audio signal from your uMic and sends the signal from your mobile phone to your hearing aids instead.


uMic and watching TV

To ensure a high quality audio in stereo sound from the TV, we recommend using the dedicated uTV 3.

In cases where an analog TV-audio output is not available or while traveling, the uMic may be used to improve the sound quality while watching TV.

For best results:

1. Try to locate the TV loudspeaker.
2. Place the uMic 20 cm (8") from the speaker.
3. Consider leaving the charger plugged in while watching TV.

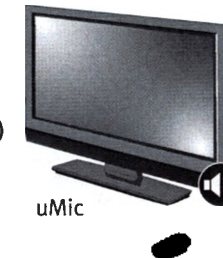uMic streaming distance

The uMic transmitter supports operation up to a range of 20 m (65 ft). Line-of-sight between your streamer and the uMic is not required. Interference in the environment by objects in the line-of-sight may reduce the operating range.

If the uMic is moved out of range of the streamer, the wireless audio transmission will be interrupted or even terminated. If the uMic is then moved back into range of the streamer, it will automatically resume streaming to your hearing aids.

It may take up to ten seconds until the streaming resumes.

 If the streaming was interrupted for more than 30 minutes, it may be necessary to switch the uMic off and on again.

Deleting all Bluetooth pairings

If you encounter issues with one of the paired devices, or if a device stops connecting to your uMic, you may consider deleting all pairings from the uMic.

1. Turn the uMic on while pushing and holding the (+) button for ten seconds. This deletes all pairings.
2. Release the (+) button when the blue indicator light turns green.

Deleting the pairing of Bluetooth devices removes the ability to connect to any previously paired Bluetooth device. You need to perform the pairing process again (see page 8) if you wish to reconnect.

 The status light will rapidly blink blue for a few seconds after pairings were deleted (automatic pairing mode).

Use of the lanyard

For transportation, you may use the detachable lanyard as illustrated.

Troubleshooting guide

Most probable cause	Remedy
I cannot hear sound from the uMic	
One or all devices are switched off	Make sure all devices (uMic, streamer, hearing aids) are turned on and fully operational
uMic is in pause mode (blue blinking status light)	Short press the streamer access button to resume streaming
The distance between the streamer and the uMic is too great	Bring devices closer to each other; 20 meters (65 ft) is the maximum distance. To avoid streaming interruptions, make sure you wear the streamer correctly and keep the uMic within the operating range
An obstacle is preventing transmission	Depending on the material, line of sight may be necessary

Most probable cause	Remedy
I hear intermittency during audio transmission	
The uMic is too far away from the streamer or shielded by an obstacle	Do not exceed the maximum distance and relocate device
The streamer searches for other paired devices for two minutes after turning on. This may cause the audio signal to drop	After two minutes, the streamer will stop searching for other devices and the intermittency will stop
The volume during streaming is not suitable (too low or too high)	
The volume levels are not set correctly	Adjust the volume using the streamer volume (+, -) buttons. The volume can also be adjusted on the uMic (+, -) buttons

Most probable cause	Remedy
There is too much noise, I cannot understand the speaker	
Your hearing aid microphones are set too loud	Reduce the hearing aid volume by pressing the streamer (M) key until you reach a comfortable level. Place the uMic closer to the speakers mouth (20 cm, 8"). Use the uMic (+) button to increase the uMic volume
The uMic rubs against the wearer's clothing	Make sure the uMic is clipped firmly. Clothing should not cover the uMic
The settings of your hearing aids are not right	Please ask your hearing healthcare professional to fine tune your hearing aids

Most probable cause	Remedy
---------------------	--------

uMic is no longer recognized by the streamer with which it was previously paired

Two minutes after no device was found, the streamer stops seeking devices to save power	Switch the streamer off. Wait until the red power indicator turns off, then turn device on again
---	--

The devices are too far away from each other	Bring devices within 1 meter range and try again
--	--

uMic can be paired to a maximum of four different Bluetooth devices. If the internal memory is full, the fifth pairing will delete the first pairing	Repeat the pairing between the devices
--	--

Pairing was never done, or it was manually deleted	Repeat the pairing between the devices
--	--

The light indicator does not turn solid blue

uMic and the streamer are not paired	Repeat the pairing between the devices
--------------------------------------	--

Most probable cause	Remedy
---------------------	--------

I see a rapidly blinking blue indicator light on the uMic

uMic is in pairing mode	Complete the pairing
-------------------------	----------------------

Indicator lights are solid blue but I cannot hear the uMic

The streamer is out of range of the hearing aids	Wear the streamer around the neck and correctly secure the neckloop
--	---

The streamer is not correctly programmed for your hearing aids	Ask your hearing healthcare professional to reconfigure your streamer and the hearing aids
--	--

The charging time is longer than usual

Device is being used during charging	Charging during use is possible, however, the charging time is increased by 30% (≥2hr)
--------------------------------------	--

The charger is not providing enough charging current	Only use chargers approved by Unitron
--	---------------------------------------

For any problems not listed in this user guide, please contact your hearing healthcare professional.

Warnings

- Please read the information on the following pages before using your Unitron uMic.

Hazard warnings

- ⚠ Keep this device out of reach of pets and children under 3 years.
- ⚠ Only use hearing aids that have been specially programmed for you by your hearing healthcare professional.
- ⚠ Changes or modifications to any of the devices that were not explicitly approved by Unitron are not permitted.
- ⚠ Opening the uMic might damage it. If problems occur which cannot be resolved by following the remedy guidelines in the troubleshooting section of this user guide, consult your hearing healthcare professional.
- ⚠ This product has an embedded, non-replaceable battery. Do not attempt to open

the product or remove the battery as this may cause injury and damage the product.

- ⚠ Dispose of electrical components in accordance with your local regulations.
- ⚠ Do not use the device in explosive areas (mines or industrial areas with danger of explosions, oxygen rich environments or areas where flammable anesthetics are handled) or where electronic equipment is prohibited.
- ⚠ External devices may only be connected if they have been tested in accordance with corresponding IECXXXXX standards. Only use accessories approved by Unitron.
- ⚠ Using your uMic cables in any way contradictory to their intended purpose (e.g., wearing the USB cable around the neck) can cause injury.
- ⚠ When operating machinery, ensure that no parts get caught in the machine.
- ⚠ The mini-USB port is to be used for the described purpose only.

- ⚠ Caution: electric shock. Do not insert plug alone into electrical outlets.

- ⚠ For safety reasons, only use chargers supplied by Unitron or stabilized chargers with a rating of 5VDC, min. 500 mA.

Information on product safety

- ⚠ Do not charge uMic from a PC or laptop USB port. This may damage your equipment.
- ⚠ Turn off the uMic when not in use.
- ⚠ Always plug in the uDirect 2 or uDirect 3 neckloop and wear it around the neck when transmitting signals to your hearing aids.
- ⚠ Protect the device connectors, plugs, charging slot and power supply from dirt and debris.
- ⚠ Protect all devices from excessive moisture (bathing or swimming areas) and heat

28

sources (radiator, TV top). Protect devices from excessive shock and vibration.

- ⚠ Clean devices using a damp cloth. Never use household cleaning products (washing powder, soap, etc.) or alcohol to clean them. Never use a microwave or other heating devices to dry any of the devices.
- ⚠ The digitally-coded, inductive transmission technology used in these devices is extremely reliable and experiences virtually no interference from other devices. It should be noted, however, that when operating the system near computer equipment, larger electronic installation or other strong electromagnetic fields, it may be necessary to be at least 60 cm (24") away from the interfering device to ensure proper operation.
- ⚠ Special medical or dental examination including radiation described below, may adversely affect the correct functioning of your uMic. Remove and keep it outside the examination room/area before undergoing

29

a medical or dental examination with X-ray (also CT scan) or any medical examinations with MRI/NMRI scans, generating magnetic fields.

Other important information

- ⚠ The uMic is a Bluetooth 2.1 compliant device. Its wireless transmission technology is not limited to but optimized for Unitron streamer devices such as the uDirect 2 or uDirect 3.
- ⚠ The use of uMic with non-Unitron streamers is not supported by Unitron.

Information and explanation of symbols

With the CE symbol, Unitron confirms that this Unitron product – including accessories – meets the requirements of the Medical Devices Directive 93/42/EEC as well as the Radio Equipment Directive 2014/53/EU on radio and telecommunications transmitters. The numbers after the CE symbol correspond to the code of certified institutions that were consulted under the above-mentioned directives.

This symbol indicates that it is important for the user to read and take into account the relevant information in this user guide.

This symbol indicates that it is important for the user to pay attention to the relevant warning notices in this user guide.

Important information for handling and product safety.

Australia and New Zealand EMC and Radiocommunications compliance label.

This symbol indicates that the products described in these user instructions adhere to the requirements for an application part of Type B of EN 60601-1. The surface of the hearing aid is specified as applicated part of Type B.

Copyright symbol

This symbol shall be accompanied by the name and the address of the authorised representative in the European Community.

This symbol shall be accompanied by the name and the address of the manufacturer (who are placing this device on the market).

Operating conditions:

This device is designed such that it functions without problems or restrictions if used as intended, unless otherwise noted in these user guides.

Transport conditions:

Temperature: -20° to +60° Celsius (-4° to +140° Fahrenheit).
Humidity: Up to 90% (non condensing).

The symbol with the crossed-out garbage bin is to make you aware that this device may not be thrown away as normal household waste. Please dispose of old or unused devices, at waste disposal sites intended for electronic waste, or give your device to your hearing healthcare professional for disposal. Proper disposal protects the environment and health.

The Bluetooth® word mark and logos are registered trademarks owned by Bluetooth SIG, Inc. and any use of such marks by Unitron is under license. Other trademarks and trade names are those of their respective owners.

uMic is a trademark of Unitron.

Compliance information

Declaration of Conformity

Hereby, Unitron Hearing declares that this product meets the requirements of the Medical Device Directive 93/42/EEC and complies with Radio Equipment Directive 2014/53/EU. The full text of the Declarations of Conformity can be obtained from the manufacturer.

Notice 1:

This device complies with Part 15 of the FCC Rules and with RSS-210 of Industry Canada. Operation is subject to the following two conditions:

- 1) this device may not cause harmful interference, and
- 2) this device must accept any interference received, including interference that may cause undesired operation.

Notice 2:

Changes or modifications made to this device not expressly approved by Unitron may void the FCC authorization to operate this device.

Notice 3:

This device has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules and ICES-003 of Industry Canada. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This device generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference

to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this device does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

Reorient or relocate the receiving antenna.

Increase the separation between the device and receiver.

Connect the device into an outlet on a circuit different from that to which the receiver is connected.

Consult the dealer or an experienced radio/TV technician for help.

Notice 4:

FCC/IC RF Radiation Exposure Statement. This equipment complies with FCC/IC radiation exposure limits set forth for an uncontrolled environment. End users must follow the specific operating instructions for satisfying RF exposure compliance. This transmitter must not be co-located or operating in conjunction with any other antenna or transmitter.

Australia: Supplier Code Number N15398

New Zealand: Supplier Code Number Z1285

Unitron Hearing, a division of
National Hearing Services Inc.
20 Beasley Drive, P.O. Box 9017,
Kitchener, ON N2G 4X1 Canada

Sonova AG
Laubisrütistrasse 28
CH-8712 Stäfa, Switzerland

For a listing of Unitron group
companies, please visit
www.unitron.com

CE 0678

Appendix 1 – 029-5905-02f_uMic_Gd_ENG

Distributor

© 2015, 2017 Unitron. All rights reserved.

unitron.com

Appendix 2 – 029-5989-02b_RemoteControl2_Gd_ENG

Remote control 2 guide

Amir Panah

unitron. Hearing matters

Thank you

Thank you for choosing remote control 2 for your Unitron hearing aids.

At Unitron, we care deeply about people with hearing loss. We work closely with hearing healthcare professionals to make advanced, purpose-driven solutions available to everyone. Because hearing matters.

The intended use of the remote control 2 is to enable easy changes of your hearing aid programs and volume.

 The remote control 2 must be configured for your hearing aids by your hearing healthcare professional before you can use it.

Your remote control 2

Hearing healthcare professional: _____

Telephone: _____

Serial number: _____

Replacement batteries: One AAA (1.5 Volt)

Warranty: _____

Date of purchase: _____

Battery size

AAA

Changing batteries

Low battery warning

Light indicator changes to yellow or red

Quick reference

- 1 On/Off switch
- 2 Indicator light
- 3 +/- buttons
- 4 Program change button
- 5 Home button
- 6 Battery compartment

On/Off


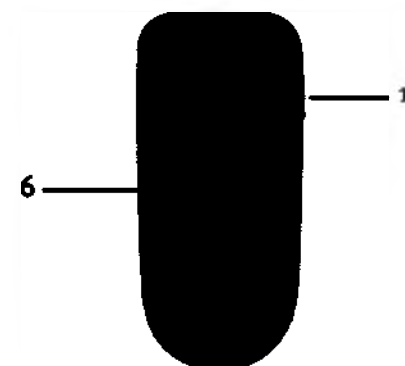
Table of contents

Your remote control 2 at a glance.....	2
Turning your remote control 2 on and off.....	4
Replacing the battery.....	5
Operating instructions.....	7
Understanding the indicator light.....	9
Adjusting your hearing aids with your remote control 2	10
Caring for your remote control 2.....	12
Cleaning your remote control 2	13
Troubleshooting guide	14
Warnings	16
Information and explanation of symbols.....	21
Compliance information.....	23

Your remote control 2 at a glance


- 1 **On-off switch** - turn your remote control on or off and lock it to avoid sending unintentional commands
- 2 **Indicator light** - will light up when you turn on the remote control or press any button
green light = battery OK
yellow light = low battery
red light = very low battery
- 3 **+/- buttons** - adjusts the volume of your hearing aids
- 4 **Program change button** - advances your hearing aid to the next listening program
- 5 **Home button** - returns your hearing aids to their start-up program and start-up settings
- 6 **Battery compartment** - allows access to change the battery

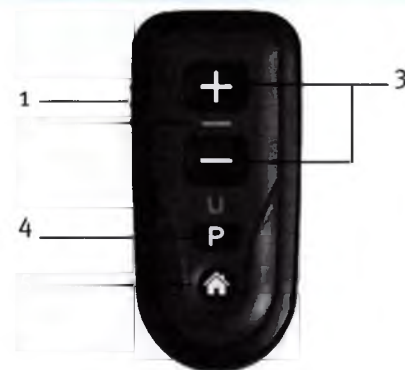Turning your remote control 2 on and off

Your remote control has an On-Off switch located on the left side of the remote.

- To turn the remote control on, slide the switch to the On or unlocked position, revealing a white dot.
- To turn the remote control off, slide the switch to the Off, hiding the white dot. Sliding the switch to the Off position will avoid sending unintentional commands to your hearing aids (while the device is in your handbag or pocket).

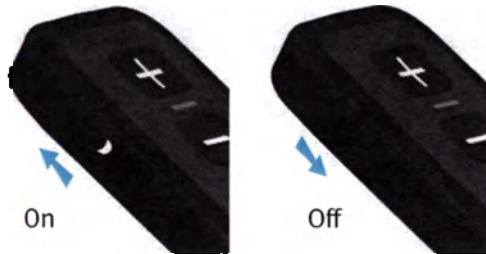4

Replacing the battery

The battery compartment is located on the back of the remote.

1. Slide the battery cover in the direction of the arrow.
2. Grasp the old battery with your thumb and index finger and remove.
3. Insert the new battery into the battery compartment in the correct +/– position, placing the battery's minus (–) end into the compartment first.
4. Slide the battery cover in the opposite direction of the arrow.

 The mini-USB port cannot be used for battery charging. It is to be used for programming purposes by your hearing healthcare professional.

5

Low battery warning

When the indicator light shines yellow or red instead of green, your remote control is warning you that the battery in your remote is low. You will see the warning every time you turn on your remote control or send a command to your hearing aids until you change the battery or the battery dies.

Caring for batteries

Always discard batteries in a safe and environmentally friendly way.

To prolong battery life, remember to turn the remote control off when not in use.

Operating instructions

Before using your remote control 2, turn on your hearing aids to ensure they are able to receive the remote commands.

Your remote is designed to be used comfortably in your hand at a maximum distance of one arm's length (1 m/40"). Hold the remote so the buttons are facing you, however, you do not have to point the remote at your hearing aids. Placing the remote directly next to or too close to your hearing aids may decrease the remote's effectiveness.

- ⚠ The mini-USB port is to be used for programming purposes only by your hearing healthcare professional.

- ⚠ Interference with medical devices:
This remote control uses low-power, digitally coded transmission to communicate to your hearing system. Although unlikely,

interference with medical devices such as pacemakers is possible. Therefore pacemaker users must not keep this remote in or near a breast pocket.

Understanding the indicator light

When turning on the remote control 2, the indicator light will switch on for two seconds, indicating that the remote is ready to use. the indicator light remains on as long as a button is pressed.

 If the indicator does not blink when a button is pressed, the remote needs to be programmed by your hearing healthcare professional.

The indicator light will be either green or red, depending on the battery status.

Green: battery OK

Red: < 10% battery left.

Replace the battery immediately.

Adjusting your hearing aids with your remote control 2

+/- buttons as volume control

Press + to increase the volume

Press – to decrease the volume

As the volume is adjusted, your hearing aids may beep. Please refer to your hearing aid guide for more information on your hearing aid volume beeps.

Program change button

Press the program change button to change to the next listening program in your hearing aids. Please refer to your hearing aid guide to determine which programs are available to you and your program beep notifications.

Home button

Press the home button to return your hearing aids and settings to their start-up program and the settings. Please refer to your hearing aid guide to determine which program is your start-up program.

Press and hold the home button for two seconds to activate your main hearing program as configured by your hearing healthcare professional.

 Your hearing aids will beep to confirm chosen selection.

 Changes are applied to both hearing aids at the same time.

Caring for your remote control 2

- Open the battery compartment and remove the battery if the remote control will not be used for long periods of time.
- Do not use your remote in the bath or shower or immerse it in water.
- If your remote does become wet, do not attempt to dry it in an oven or microwave. Do not adjust any controls. Open the battery compartment and allow your remote control to dry naturally for 24 hours.
- Protect your remote from excessive heat (hair dryer, vehicle glove box or dashboard).
- Do not drop your remote or knock it against hard surfaces.
- If the hearing aids do not respond to the device because of an unusual field disturbance, move away from the disturbing field.

Cleaning your remote control 2

Use a damp cloth to clean the outside of the remote. Do not use household cleaning products (washing powder, soap, etc.) or alcohol to clean the device. Never use a microwave or other heating device to dry your remote.

Troubleshooting guide

You can refer to the troubleshooting table below if you are experiencing problems with your remote control. For problems with your hearing aids, please consult your hearing aid guide.

Cause	Possible remedy
Indicator light not working	
Not turned on	Turn on
Low/dead battery	Replace battery
Poor battery contact	Consult your hearing healthcare professional
Battery upside down	Insert battery plus (+) side up
Indicator light is yellow or red	
Low battery	Replace battery
Hearing aids not responding to remote control commands	
Holding the remote too far away	Move the remote a maximum distance of one arm's length

For any problems not listed in the guide, contact your hearing healthcare professional.

Warnings

 Please read the information on the following pages before using your Unitron remote control.

Hazard warnings

-  Keep this device out of reach of pets and children under 3 years.
-  Only use hearing aids that have been specifically programmed for you by your hearing healthcare professional.
-  Changes or modifications to any of the devices that were not explicitly approved by Unitron are not permitted.
-  This remote control uses low-power, digitally coded transmission to communicate to your hearing system. For persons with active implantable medical devices, (i.e. pacemakers, defibrillators, etc.) caution should be exercised when using this device.
-  Dispose of electrical components in accordance with your local regulations.

- ⚠ Do not use the device in explosive areas (mines or industrial areas with danger of explosions, oxygen rich environments or areas where flammable anesthetics are handled) or where electronic equipment is prohibited.
- ⚠ External devices may only be connected if they have been tested in accordance with corresponding IECXXXXX standards. Only use accessories approved by Unitron.
- ⚠ When operating machinery, ensure that no parts get caught in the machine.
- ⚠ The mini-USB port is to be used for the described purpose only.
- ⚠ Using the device, its connectors or cables in any way contradictory to their intended purpose as described in this user guide (e.g., wearing the USB cable around the neck) can cause injury.
- ⚠ Only use accessories approved by Unitron to avoid any electrical shock.

Battery warnings

- ⚠ Never leave your remote control or batteries where small children and pets can reach them.
- ⚠ Remove the battery if you are not using the remote control for a long period of time.
- ⚠ Never put your remote control or batteries in your mouth. If a battery is swallowed, call a physician immediately.
- ⚠ Always use a high quality, long-life alkaline battery as a replacement and dispose of the old one in an environmentally-friendly way. Do not use rechargeable batteries.

Precautions

- ⚠ Your remote control should only be used as directed by your hearing healthcare professional.
- ⚠ If your remote control fails to operate after you have correctly inserted a new battery, contact your hearing healthcare professional.

- ⚠ Special medical or dental examination including radiation described below, may adversely affect the correct functioning of your remote. Remove and keep it outside the examination room/area before undergoing a medical or dental examination with X-ray (also CT scan) or any medical examinations with MRI/NMRI scans, generating magnetic fields.
- ⚠ Always bring both your remote control and your hearing aids for service questions.
- ⚠ Your remote control and your hearing aids are given a unique communication code. Because of this code, your remote control will not affect hearing aids worn by others.
- ⚠ Your remote control uses the most modern components to provide an extremely reliable experience with virtually no interference from other devices. However, when operating the remote control near a computer terminal, communication device

such as a digital cell phone, or other strong electromagnetic field, it may be necessary to be at least 24 inches (60 centimeters) away to ensure proper operation.

- ⚠ If your hearing aids do not respond to the remote control because of an unusual field disturbance, move away from the disturbing field.
- ⚠ Keep the device at least 10 cm away from any kind of magnets.
- ⚠ High-powered electronic equipment, larger electronic installations and metallic structures may impair and significantly reduce the operating range.

Labeling

The serial number and year of manufacture are located inside the battery compartment.

Please fill out the contact information sticker provided with your new remote control and affix it to the inside of the battery door.

Information and explanation of symbols

With the CE symbol, Unitron confirms that this Unitron product – including accessories – meets the requirements of the Medical Devices Directive 93/42/EEC as well as the Radio Equipment Directive 2014/53/EU on radio and telecommunications transmitters.

The numbers after the CE symbol correspond to the code of certified institutions that were consulted under the above-mentioned directives.

This symbol indicates that it is important for the user to read and take into account the relevant information in this user guide.

This symbol indicates that it is important for the user to pay attention to the relevant warning notices in this user guide.

Important information for handling and effective use of the product.

Australia and New Zealand EMC and Radiocommunications compliance label.

This symbol indicates that the products described in these user instructions adhere to the requirements for an application part of Type B of EN 60601-1. The surface of the hearing aid is specified as applicated part of Type B.

Copyright symbol

This symbol shall be accompanied by the name and the address of the authorised representative in the European Community.

This symbol shall be accompanied by the name and the address of the manufacturer (who are placing this device on the market).

Operating conditions:

This device is designed such that it functions without problems or restrictions if used as intended, unless otherwise noted in these user guides.

-10° – 55° Celsius and relative humidity of <95% (non condensing).

Transport conditions:

Temperature: -20° to +60° Celsius (-4° to +140° Fahrenheit).

Humidity: Up to 90% (non condensing).

The air pressure between 500 and 1100 hPa is appropriate.

The symbol with the crossed-out garbage bin is to make you aware that this device may not be thrown away as normal household waste. Please dispose of old or unused devices, at waste disposal sites intended for electronic waste, or give your device to your hearing healthcare professional for disposal. Proper disposal protects the environment and health.

Compliance information

Declaration of Conformity

Hereby, Unitron Hearing declares that this product meets the requirements of the Medical Device Directive 93/42/EEC and complies with Radio Equipment Directive 2014/53/EU. The full text of the Declarations of Conformity can be obtained from the manufacturer.

Notice 1:

This device complies with Part 15 of the FCC Rules and with RSS-210 of Industry Canada. Operation is subject to the following two conditions:

- 1) this device may not cause harmful interference, and
- 2) this device must accept any interference received, including interference that may cause undesired operation.

Notice 2:

Changes or modifications made to this device not expressly approved by Unitron may void the FCC authorization to operate this device.

Notice 3:

This device has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules and ICES-003 of Industry Canada. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This device generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this device does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the device and receiver.
- Connect the device into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

Australia: Supplier Code Number N15398

New Zealand: Supplier Code Number Z1285

Appendix 2 – 029-5989-02b_RemoteControl2_Gd_ENG

Appendix 2 – 029-5989-02b_RemoteControl2_Gd_ENG

Unitron Hearing, a division of
National Hearing Services Inc.
20 Beasley Drive, P.O. Box 9017,
Kitchener, ON N2G 4X1 Canada

Sonova AG
Laubisrütistrasse 28
CH-8712 Stäfa, Switzerland

For a listing of Unitron group
companies, please visit
www.unitron.com

Appendix 2 – 029-5989-02b_RemoteControl2_Gd_ENG

Distributor

8/17/2016 00:00:00 Unitron. All rights reserved.

unitron.com

Appendix 3 – 029-5990-02b_uDirect3_Gd_ENG

uDirect3 guide

Amir Panah

unitron. Hearing matters

Thank you

Thank you for choosing the Unitron uDirect™ 3. At Unitron, we care deeply about people with hearing loss. We work closely with hearing healthcare professionals to make advanced, purpose-driven solutions available to everyone. Because hearing matters.

The intended use of uDirect 3 is to wirelessly stream audio from different Bluetooth® devices to your hearing aids. For example, mobile phones, MP3 players etc.

You can also stream other audio devices directly via a plug. For example, computers, Roger and FM receivers.

uDirect 3 can also be used as a remote control to adapt your hearing aids to different listening situations.

 Please note that your hearing healthcare professional must configure your uDirect 3 before it can be used as a remote control for your hearing aids.

Your uDirect 3

Hearing healthcare professional: _____

Telephone: _____

Serial number of uDirect 3: _____

Warranty: _____

Date of purchase: _____

The uControl™ app can be downloaded at the App Store and Google play. Search for Unitron and select uControl.

Quick reference


Charging uDirect 3


Wearing uDirect 3

Bluetooth pairing

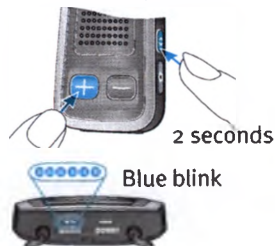On/Off

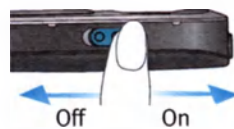Power indicator

Solid red - charging
Solid green - fully charged
Red blink - low battery
No light - battery drained

Audio indicator

Solid orange - audio jack or FM streaming
Solid violet - demo sound
Solid blue - Bluetooth streaming
Blue short blink - 1 device connected
Blue double blink - 2 devices connected
Blue very fast blink - pairing

Table of contents

Your uDirect 3 at a glance	2
Setting up the power supply	4
Charging the uDirect 3	5
Turning the uDirect 3 on and off	6
Summary of indicator lights	7
Wearing the uDirect 3 correctly	10
Operating distance	11
Button lock function	11
Testing sound via the uDirect 3	12
Access button functionality	14
Spoken notifications	15
uControl app	16
Using uDirect 3 as a remote control for your hearing aid	17
Using uDirect 3 with a mobile phone	20
Advanced phone options	31
Using the external microphone	35
Using uDirect with uTV 3	36
Using uDirect 3 with uMic	40
Using uDirect 3 with an audio device with Bluetooth	42
Using uDirect 3 with an audio device with a headphone jack ..	45
Setting up the FM receiver	46
Managing multiple audio inputs to the uDirect 3	49
Resetting your uDirect 3	52
Deleting the pairing of Bluetooth devices	53
Troubleshooting guide	54
Warnings	64
Information and explanation of symbols	71
Compliance information	73

Your uDirect 3 at a glance

Buttons

- 1 Access (A)
- 2 Volume (+/-)
- 3 Home (house icon)
- 4 Connect (hand icon)
- 5 On/Off switch (I/O)

Inputs

- 6 Microphone openings
- 7 Audio input (3.5 mm jack)
- 8 FM input (Europlug)
- 9 Charging input (mini-USB)


Indicators

- 10 Audio indicator
- 11 Power indicator (power)

Neckloop


- 12 Neckloop sockets
- 13 Neckloop (with antenna)
- 14 Neckloop plugs

*Can be configured by your hearing healthcare professional.

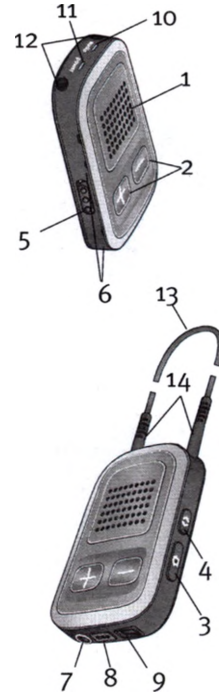Accessories

- 15, 16 Universal power supply with mini-USB plug
- 17-20 US, EU, UK, AUS adapter

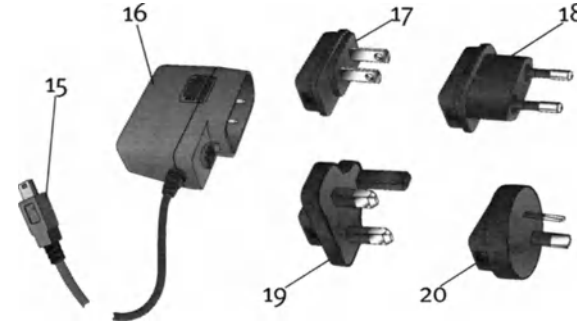Other optional accessories

- External lapel microphone 'MC1'
- Short neckloop (51 cm)
- Unitron uTV 3 (Bluetooth transmitter for TV use)
- Audio cable, 1.5 m / 5 ft.
- USB cable, 3 m / 10 ft.
- (Included country adaptors may vary by country.)

Setting up the power supply

1. Select the adapter that corresponds to your country.

2. Insert the round corner of the adapter to the round corner of the universal power supply until it is completely inserted.

3. Click the tail of the adapter in to lock it into place. Make sure the adapter is securely locked.

To remove the power adapter:

1. Pull back the tab on the universal power supply and hold.
2. Gently pull the adapter upwards to remove it from receptacle.
3. Then start over with step one (above).

Charging the uDirect 3

Your uDirect 3 has a non-removable rechargeable battery.

To charge it:

1. Plug the mini-USB plug into the mini-USB port on the uDirect 3.

2. Plug the power supply into a power outlet. During charging, the power indicator is solid red, this typically takes 90 minutes.

When the battery is fully charged, the power indicator is solid green.


It is safe to leave the uDirect 3 connected to the charger overnight.

 When charging your uDirect 3 for the first time, let it charge for at least three hours, even if the battery indicator is not lit or turns green earlier.

 uTV™ 3 can also be used as a charging station for uDirect 3.

Turning the uDirect 3 on and off

Slide the power switch (I/O) as shown to turn the uDirect 3 on or off.

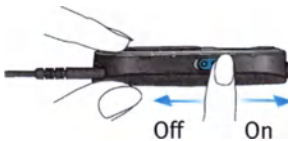Switch the uDirect 3 off when not in use or if in storage or during transportation.

For operational use switch the uDirect 3 on.

Summary of indicator lights

The indicator lights provide helpful information as follows:

Power indicator

Charging	Solid red	
Charged and full battery	Solid green	
Switching on	Green 2 seconds	
Switching off	Blink red, then solid 2 seconds*	
Battery above 10%, normal operation	Slow blink green **	
<10% battery, charging required (at least 2 hours)	Short blink red**	
Battery dead, charging required (at least 2 hours)	No light after switching on	
Enable/disable lock state	Short blink red (lock) or short blink green (unlocked)	

- * When switching the uDirect 3 off, the power indicator may blink red first while the uDirect 3 is properly disconnecting all Bluetooth devices.

** If the neckloop is not connected to the uDirect 3, the battery level is only shown during startup, charging and when a button is pressed.

 After unplugging the neckloop the uDirect 3 enters power down mode after two minutes and the power indicator goes off. Only the remote control functions can be used while the neckloop is unplugged.

Battery state

The color of the power indicator informs you about the uDirect 3 battery state.

The green power indicator indicates a battery level >10%.

The red power indicator indicates a battery level <10%. Less than one and a half hours streaming remain and the battery should be charged for at least two hours.

Audio indicator

The neckloop must be plugged in for any streaming.

Plug or FM streaming	Solid orange	
Bluetooth phone or music streaming	Solid blue	
One Bluetooth device connected	Short blink blue	
Two Bluetooth devices connected	Double-blink blue	
Bluetooth pairing	Very fast blue blinking	
Demo sound playing	Solid violet	

 The Bluetooth pairing indication will be visible even if the neckloop is not plugged in.

 During cable or FM streaming, the audio indicator will remain orange regardless of the Bluetooth connection state.

Wearing the uDirect 3 correctly

Before any phone, music or TV streaming can take place, the uDirect 3 must be switched on and worn correctly as shown below:

1. Put the neckloop around your neck (easily done by unplugging one neckloop plug).
2. Make sure both plugs are well connected in the neckloop sockets.
3. Switch the uDirect 3 on.

 For streaming, the neckloop must always be plugged in. When unplugged, the uDirect 3 will go into sleep mode after two minutes. The power and audio indicators will stop blinking and all active Bluetooth connections are terminated to save power. Plug the neckloop in again to activate the streaming.

Operating distance

If uDirect 3 is used as a handheld remote control (without neckloop):

- Maximum distance 110 cm/43"
- Hold uDirect 3 as shown

Button lock function

The access (A) button on your uDirect 3 can be locked to avoid inadvertent presses. To lock the button, press and hold the connect (↔) and Home (🏠) buttons simultaneously. The power indicator light will flash red three times to indicate that the buttons have been locked.

Press and hold the connect (↔) and Home (🏠) buttons simultaneously again to unlock the buttons.

You will still be able to answer an incoming telephone call using the access (A) button while it is locked.

-  The power indicator is lit as long as you press a uDirect 3 button. Green confirms a good uDirect 3 battery state, red indicates that you need to charge the uDirect 3 battery.
-  Two minutes after unplugging the neckloop both indicator lights will switch off and the uDirect 3 disconnects Bluetooth to save power.

Testing sound via the uDirect 3

1. Wear the uDirect 3 correctly.
2. Press and hold the access (A) button while switching on the device. Release the access button. The audio indicator light will turn violet.

The hearing aids will switch to the streaming mode and you will hear the uDirect 3 sound demo.

Use the volume (+/-) buttons to adjust the volume.

The sound demo will stop after three minutes or when the access button is pressed.

Adjusting the volume during the sound demo

will also set the volume level of the spoken uDirect 3 messages.

If you cannot hear the sound demo:

- Make sure your hearing aids are inserted correctly.
- Make sure the uDirect 3 neckloop is properly plugged in on both ends and worn around your neck.
- Make sure the uDirect 3 is charged and switched on.
- Make sure you are not on a phone call connected through the uDirect 3.

Access button functionality

The access button has multiple functions:

- **Phone:** Manages telephone functions such as accepting and ending phone calls. See the section “Using uDirect 3 with a smartphone or cell phone” for more information.
- **Streaming:** It pauses and resumes the streaming from an audio device.
- **Remote control:** It can be used to select the hearing aid programs.

Spoken notifications

Spoken notifications inform you about the different operation modes and status of the uDirect 3 using text-to-speech, for example low battery.

Your hearing healthcare professional can configure the language setting as well as which notifications are spoken.

- 📖 **Caller identification is independent of the spoken notifications level.**
- 📖 **Spoken notifications must be enabled by your hearing healthcare professional.**

uControl app

The uControl is a smartphone app which offers quick and easy volume and program adjustments.

It also provides:

- individual left and right volume control
- additional program controls
- easy Bluetooth pairing


Using uDirect 3 as a remote control for your hearing aid

 Before the remote control functionality can be used, your hearing healthcare professional must enable the remote control.

There are two ways to use the uDirect 3 as a remote control:

- The neckloop is plugged in with both ends, and the uDirect 3 must be worn around your neck. The remote control function is available as long as you are not streaming audio.
- The neckloop is completely un-plugged and uDirect 3 is operated at arm's length. No streaming is possible in this configuration.

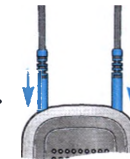Make sure the uDirect 3 is switched on.

Volume change

Press the (+/−) volume button on the uDirect 3 to adjust the volume of your hearing aids:

- Press (+) to increase loudness
- Press (−) to decrease loudness

Program change

Use the access (A) button to change hearing aid listening programs.

Changing hearing aid programs is only possible if no streaming is in progress.

Press and hold the home (🏠) button (> two seconds) to select your favorite hearing program configured by your hearing healthcare professional.

Push the home (🏠) button to return to the startup hearing aid program and volume levels of your hearing aids.

More remote control information

- Your hearing aids may beep to confirm that your selection was applied.
- Both hearing aids are controlled at the same time.
- The available hearing aid programs and volume levels are set by your hearing healthcare professional.

Using uDirect 3 with a mobile phone

If you have a mobile phone equipped with Bluetooth functionality (e.g. iPhone or android), uDirect 3 may be used to stream phone conversations wirelessly to your hearing aids allowing you to hear the callers voice directly in your hearing aids. uDirect 3 picks up your voice through its own microphone and streams it to your mobile phone.

During phone calls, the access (A) button provides a range of features that help you manage your calls.

To benefit from the pairing of the uDirect 3 and your mobile phone, we recommend wearing the uDirect 3 around the neck. During calls, your voice can only be picked up if the uDirect 3 is not covered by clothes.

If you wear the uDirect 3 under clothing during phone calls, the external Unitron microphone “MC1” is needed (see the section in this guide on using the external microphone).

Pairing your mobile phone

1. Enable Bluetooth on your mobile phone settings menu.
2. Ensure that the uDirect 3 is turned on and the neckloop is plugged in and worn correctly. Press the connect (↔) and the volume up (+) buttons at the same time for two seconds. The audio indicator will blink blue to indicate pairing mode on the uDirect 3.

3. Enable Bluetooth pairing mode on your mobile phone so that a new device can be added. See manufacturer instructions for details. If prompted, enter “0000” for pairing code.
4. If prompted, accept request to grant uDirect 3 access to your phone book. It is used for voice calling and caller ID features.
5. Confirmation of successful connection is indicated on the uDirect 3 when the audio indicator stops blinking blue. Additionally a headset symbol (🎧) may appear on the main screen of your mobile phone.

Once pairing is done, your mobile phone should pair automatically with the uDirect 3.

-  Please contact your hearing healthcare professional if you can not complete the pairing procedure with the steps described.
-  The pairing process does not require the neckloop to be plugged in. For streaming and Bluetooth connection, the neckloop is mandatory.

-  To receive phone calls your phone must support either the Bluetooth “Head Set” (HSP) or “Hands Free” (HFP) profile. When pairing to your phone, it may ask you which Bluetooth service you want to enable. Select “Headset” or, if available, “Stereo Headset.”
-  To access caller ID features, your phone must support the Phone Book Access Profile (PBAP)
-  When pairing to a cellular phone, we recommend disabling the phone key tones and SMS alerts.
-  To listen to stereo music from your phone, it must support the A2DP/“Stereo Headset” profile.
-  A maximum of eight devices can be paired to the uDirect 3. A ninth pairing will overwrite the pairing of the device that has not been used in the longest amount of time.
-  The uDirect 3 can support (connect) two phones at the same time but only one call at a time is possible.


Preparing for phone calls

The Bluetooth range (maximum distance between uDirect 3 and the phone) is 5-10 m (15-30 ft).

Outside this range, the connection will be lost or you will encounter intermittency during a call.

When engaged in any phone call, wear the uDirect 3 normally as shown below (figure a).

- In very noisy environments, you may move the uDirect 3 closer to your mouth (figure b).
- Avoid covering the small microphone openings on the side with your fingers or clothes (figure c) and do not speak from the front (figure d).

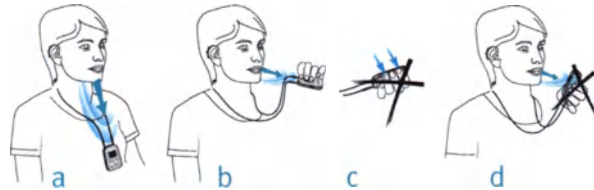Some phones are configured to play key-press tones, SMS announcements and confirmation beeps to the uDirect 3.

If this is set you will hear these tones interrupting any active audio streaming through your hearing aids. Your hearing aids may briefly switch to streaming mode which may be unwanted.

Please refer to your phone user guide for information about disabling these sounds.

- ☞ Phone calls via Bluetooth will always take priority over other audio sources on your uDirect 3.
- ☞ It is not possible to change hearing aid programs using uDirect 3 during a phone call.
- ☞ You will hear the caller's voice through both hearing aids.

Making phone calls

Use the phone keypad to enter the number and press the phone dialing button.

When connected, you will hear the dialing tone through your hearing aids after approximately three seconds.

Using cordless phones (DECT) or older phones with headset profile only

If your phone only supports the headset profile (HSP), the connection must be initiated manually for outgoing calls. This often applies to cordless home phones (DECT standard).

 For these phones, the audio indicator might not blink blue even if the phone is paired.

1. Dial the number you want to dial on your phone, or select it from the phone book.
2. Briefly press the connect () button, then wait five seconds.
3. The dialing process will automatically start and the audio indicator will turn solid blue. You may hear a confirmation beep through your hearing aids.

Incoming calls

The uDirect 3 must be switched on, paired to the phone and worn correctly.

When a call is received on your phone, the uDirect 3 audio indicator will show a solid blue light, and you will hear the ringtone signal through your hearing aids.

Accepting incoming calls

To accept an incoming call, briefly press the access () button on your uDirect 3 once you hear the ringtone through your hearing aids. It will take about five seconds before you hear the ringtone through your hearing aids.

Speak into the uDirect 3, not into the phone, when the call is in progress.


 For your safety, the hearing aid microphones will remain on at reduced levels during phone calls and audio streaming.

 Only one phone call at a time is possible. A second ringing phone is ignored while another call is already connected. You must disconnect the first call before answering the second one.

Ending a call

To end a phone call via your uDirect 3, briefly press the access () button. The call will terminate and your hearing aids will return to the previously-active program or audio source (such as the audio jack input).

You can also end a call by operating your phone instead of your uDirect 3.

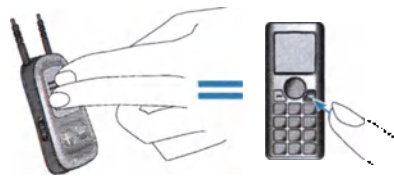Different ringtones of the uDirect 3

The phone ringtone depends on the uDirect 3 configuration and your phone capabilities.

The ringtone is assigned automatically during the initial phone pairing process and cannot be configured.

When pairing more than one phone, the second and third paired phone will have a different ringtone. A total of three different ringtones are automatically assigned.

The uDirect 3 will also announce the caller name as shown on the phone display (caller identification using text-to-speech technology). The name is taken from your phone's directory.

 Please note that this feature must be supported by your phone. Please consult the phone user guide or the phone supplier for more information.

 The default language setting of a new uDirect 3 is US English. Your hearing healthcare professional may configure the appropriate language setting for proper pronunciation of the names. Your hearing healthcare professional may switch the caller identification off upon your request.

Rejecting a call

If you prefer not to answer an incoming call, press and hold the access (A) button on your uDirect 3 for two seconds until the ringtone stops. Using your phone to reject the call will have the same effect.

Pressing the home (H) button will reject the call and the hearing aids will switch to the startup program.

Advanced phone options

Your hearing healthcare professional may configure the following additional features on uDirect 3. They will be assigned to a long press (two seconds) of the access (A) or the connect (C) button:

 Your mobile phone must support these features. Please consult the phone user manual for details.

Redial last number

A long press of the access (A) button will redial the last dialed number.

Voice dial

A long press of the connect (C) button will start the voice dial which allows you to speak the person's name instead of typing the number.

 When using two phones at the same time, only the last paired phone can be used for voice dial.

Hold call

While in a phone conversation, a long press of the connect (↔) button will put the call on hold. The caller stays on the line but there is no audio transfer and the hearing aids switch to the previous hearing program. Use this mode to have a side conversation with a person next to you without transmitting the conversation to the person on the phone.

Transfer call

While in a phone conversation, a long press of the access (A) button will transfer the call to the phone. This allows you to hand your phone to another person. Your hearing aids will switch back to the previous hearing program.

After holding or transferring a call, another long press of the access (A) button transfers the call back to the uDirect 3 and the call is resumed through your hearing aids.

 As long as the call is on hold or transferred, you may use the access (A) button to change your hearing aid programs again.

Caller identification, spoken caller names

The uDirect 3 may be able to tell you the name of the caller.

If your phone supports Phone Book Access (PBA) or Phone Book Access Profile (PBAP), the uDirect 3 will try to access information about the caller every time the phone rings.

As the phone rings, the displayed information is transferred to the uDirect 3 and spoken through your hearing aids as follows:

- The stored name in your phone's phonebook. In the example, you will hear "Mike Office"
- The phone number, if the name is not known

To make sure the pronunciation is correct your hearing healthcare professional can configure the language settings of your uDirect 3.

The phonebook is not downloaded to the uDirect 3. This way, uDirect 3 can also support the use of different phones.

Please consult your phone's user guide or the phone dealer to check whether this feature is supported by your phone.

 Some hearing aid models do not support caller identification.

 Caller identification must be enabled by your hearing healthcare professional.

Using the external microphone


Use the optional external lapel microphone “MC1” if you wear the uDirect 3 underneath clothing during a phone call.

1. Plug the 3.5mm microphone plug into the audio input of the uDirect 3, disabling the internal uDirect 3 microphones.

2. Use the clip to place the omnidirectional microphone as close to your mouth as possible.

3. Put the neckloop around your neck and switch the uDirect 3 on.


Jackets, etc., will not negatively influence the wireless link to your hearing aids.

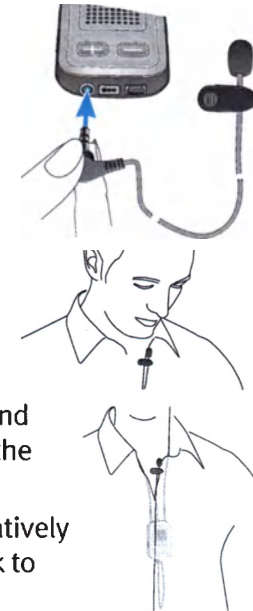Using uDirect with uTV 3

If you have a uTV~ 3, you can:

- Stream audio from your TV up to 30 m (90 ft) to your hearing aids
- Recharge your uDirect 3 in the charging cradle

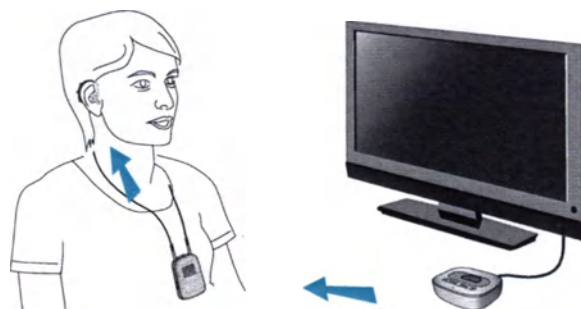Charging your uDirect 3 on the charging cradle

1. Ensure the uTV 3 has the correct charging cradle for the uDirect 3 – see uTV 3 manual for details.
2. Ensure the uTV 3 is turned on and place the uDirect 3 into the charging cradle.

3. Allow two hours for a full charge.

Pairing uTV 3 to uDirect 3

1. Switch on the uDirect 3 first. Press the connect (↔) and volume up (+) buttons at the same time for two seconds.
2. Switch on the uTV 3 and then press the pair button on the back of the uTV 3. The main indicator light on the uTV 3 will start rapidly blinking blue. Make sure uTV 3 and the uDirect 3 are within 1 m (3 ft) range.
3. After successful pairing, the blue audio indicator on the uDirect 3 and uTV 3 will stop blinking.

Playing back audio from your TV

1. Begin playback on your TV.
2. Adjust the volume level on the uTV 3 to a medium to high level.
3. Make minor volume adjustments with the uDirect 3 volume buttons (+/-).

Pausing playback

If you want to temporarily stop streaming to hear a conversation, press the (A) button to pause the streaming. Press the (A) button again to resume streaming.

Stopping playback

To stop playback on your device and switch back to regular hearing programs, press and hold the access (A) button for two seconds.

 Changing hearing aid programs is not possible while streaming.

 Please contact your hearing healthcare professional if you can not complete the pairing procedure with the steps described.

 uDirect 3 is capable of streaming up to 20 hours with a fully charged battery.

Adjusting the volume

Press the uDirect 3 volume (+/-) buttons to adjust:

- Press (+) to increase the volume
- Press (-) to decrease the volume

 Adjusting the volume on the uDirect 3 will have the same effect as adjusting the volume of the hearing aids and may also amplify surrounding sounds.

Changing the program


Changing hearing aid programs is only possible if no streaming is in progress.

Press the access (A) button to change hearing aid programs.

Press the home (H) button to return to the startup program and volume level of your hearing aids.

Using uDirect 3 with uMic

If you have a uMic, it works with uDirect 3 to stream audio from a person or device to your hearing aids.

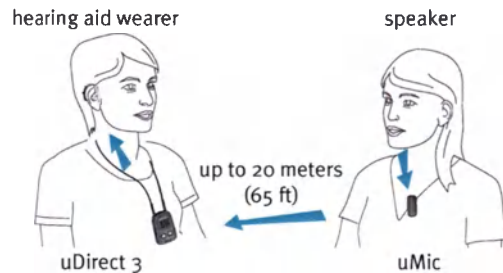Pairing uMic to uDirect 3

1. Starting with the uMic in the off position, press and hold the volume up (+) button and turn the uMic on. Continue to hold the volume up (+) button until the blue audio indicator starts to flash.
2. Switch on uDirect 3 first. Press the connect (↔) and volume up (+) buttons at the same time for two seconds.
3. After successful pairing, the blue audio

indicator on the uDirect 3 and uMic will stop blinking.

Playing back audio from your uMic

1. Clip the uMic to a person, or position it near a sound source such as a speaker or TV.
2. Adjust the volume level on the uMic to a medium to high level.
3. Make minor volume adjustments with the uDirect 3 volume buttons (+/-).

Pausing playback

If you want to temporarily stop streaming, press the (A) button to pause the streaming. Press the (A) button again to resume streaming.

Stopping playback

To stop playback on your device and switch back to regular hearing programs, press and hold the access (A) button for two seconds.

- Changing hearing aid programs is not possible while streaming.
- Please contact your hearing healthcare professional if you can not complete the pairing procedure with the steps described.

 uDirect 3 is capable of streaming up to 20 hours with a fully charged battery.

Using uDirect 3 with an audio device with Bluetooth

If you have an audio source equipped with Bluetooth functionality (e.g., MP3 player), you can stream audio wirelessly to your hearing aids.

Pairing your Bluetooth audio device:

1. uDirect 3 must be switched on and worn correctly.
2. Press the connect () and volume up () buttons at the same time for two seconds until the audio indicator starts to rapidly blink blue.
3. Start the pairing process on the Bluetooth device within one minute according to the device's user guide and use code "0000" if prompted.

4. Switch on the other Bluetooth device.
5. For some Bluetooth devices it may be necessary to press "play" on the device.
6. If the devices were paired and connected correctly, streaming will start automatically and the audio indicator will turn solid blue.

Playing back audio from your device

1. Begin playback on your device.
2. Adjust the volume level on the device to a medium to high level.
3. Make minor volume adjustments with the uDirect 3 volume buttons ( .

Pausing playback

If you want to temporarily stop streaming to hear a conversation, press the () button to pause the streaming. Press the () button again to resume streaming.

Stopping playback

To stop Bluetooth streaming completely, switch the other Bluetooth device off or disable its Bluetooth feature.

-  **Note:** Please refer to the device's user guide or dealer support if you cannot complete the pairing procedure with the generic steps described.
-  **Note:** uDirect 3 streaming time is more than 20 hours with a fully charged battery.
-  **Note:** Changing hearing programs is not possible while streaming.

Using uDirect 3 with an audio device with a headphone jack

To benefit from the streamer and remote control functions, we recommend wearing the uDirect 3 around the neck.

Streaming music via a cable

The uDirect 3 must be switched on and worn correctly around your neck.

1. Insert the audio cable into the headphones socket on your audio player and start the player.
2. Insert the other end of the audio cable into the audio input on your uDirect 3. The audio streaming is automatically activated.

Briefly press the access (A) button to stop and restart streaming.

Setting up the FM receiver

uDirect 3 can be used together with FM equipment.

1. Place the FM transmitter near the sound source and turn it on.
2. Turn the FM receiver on if required.
3. uDirect 3 must be switched on and worn correctly.
4. Attach your FM receiver to the uDirect 3 via the FM input as shown in the picture.

When FM is detected, you will hear the FM notification tone for streaming through your hearing aids. The FM notification tone is similar to a phone ringtone.

Options for handling the FM signal during the FM notification tone:

- **Accept:** Briefly press the access (A) button.
- **Put on hold:** Press the access (A) button for two seconds.

- **Ignore the FM notification tone:** After 20 seconds the FM will be put on hold automatically.

Briefly press the access (A) button to pause the FM signal during streaming.

The FM input that is paused is automatically kept in the background and can be restarted as long as an FM signal is present by briefly pressing the access (A) button.

As long as an FM signal is streaming the hearing aid programs cannot be changed via the access (A) button.

Note: After 20 seconds the input will be put on hold automatically if you do not accept it.

If uDirect 3 does not receive a signal through the Europlug for three minutes, an alert will be heard. This alert will be repeated every 15 minutes until you stop the Europlug streaming.

The uDirect 3 will play a confirmation sound each time the streaming via FM is activated and the audio indicator will turn orange.

The uDirect 3 will send a short alert signal during FM streaming if there is no FM signal for 45 seconds, and then every 15 minutes thereafter.

If there is no FM signal for more than 15 minutes, the notification tone restarts as soon as the FM signal is detected again.

To deactivate FM completely, unplug the FM receiver from the uDirect 3.

Managing multiple audio inputs to the uDirect 3

When an audio source is connected via the uDirect 3, your hearing aids will automatically switch to special programs.

It is possible to have multiple audio sources connected and in operation with the uDirect 3:

- Two Bluetooth devices (e.g., phone and uTV 3)
- One audio input (3.5 mm audio plug)
- One direct audio input (Europlug)

Spoken status indicators and the audio indicator light inform you about the connected inputs.

An audio source is considered as “connected” to the uDirect 3 if:

- The audio plug is inserted (with or without audio signal)

- An FM receiver is inserted and receiving an FM signal
 - A Bluetooth music source is streaming music or paused
 - A Bluetooth phone is paired and within reach
-  A phone call always has priority. You will hear the ringing indication even when listening to other audio sources via cable, Bluetooth or FM. If you reject or terminate the call, audio streaming will resume.
-  It is not possible to change hearing aid programs using the access (A) button as long as any audio source is connected to the uDirect 3.
-  If you want to listen to music from your connected mobile phone, you must start the music via your phone initially; it cannot be started from the uDirect 3.

Selecting from multiple inputs using the uDirect 3

With only one audio source (e.g., from a Bluetooth transmitter), each press of the access (A) button will pause or resume the audio source.

Multiple sound sources can be accessed sequentially by pressing the access (A) button. After each press, the uDirect 3 selects the next connected source in a fixed sequence.

If a source is not connected, the uDirect 3 will skip over it and proceed to the next source.

Every connected audio source has a fixed place in the sequence, as shown in the diagram.

-  Pressing the home (H) button will always select the startup programs of your hearing aids, regardless of the audio source sequence.

Resetting your uDirect 3

If, for some reason, the uDirect 3 becomes unresponsive:

1. Slide the uDirect 3 power switch to the on position.
2. Press and hold the (↔), (⏻) and (−) button simultaneously for two seconds.
3. Slide the uDirect 3 power switch to the off position.
4. Wait for five seconds.

The uDirect 3 is now reset and will resume normal operation when switching on again.

 After resetting the Bluetooth pairings, the configuration status and the grouping to the hearing aids are kept.

Deleting the pairing of Bluetooth devices

If you encounter issues with one of the paired devices, or if a device stops connecting to your uDirect 3, you may consider deleting the pairing of devices to the uDirect 3.

1. Turn on your uDirect 3.
2. Press and hold the (↔) and (+) buttons for ten seconds until the blue rapid blinking audio indicator stops.
3. Also delete the uDirect 3 pairing entries on your mobile phone before any re-pairing.
4. Repeat the pairing procedure.

 Deleting the pairing of Bluetooth devices removes the ability to connect to these devices. You will need to perform the pairing process again if you wish to reconnect.

Troubleshooting guide

Cause	Possible remedy
During Bluetooth pairing my other device is prompting for a code, what do I have to enter?	
The pairing code is required for the initial coupling of uDirect 3 with your Bluetooth device (pairing)	The uDirect 3 pairing code is "0000" (four zeros)

I hear intermittency during audio streaming via uDirect 3

The Bluetooth device is too far away, or worn in a back pocket (body shielding)	Do not to exceed the maximum distance 5-10 m (15-30 ft) and keep Bluetooth devices in front of you
uDirect 3 searches for known devices for two minutes after switching on	After two minutes, uDirect 3 will stop searching and the intermittency will stop
A phone connected to other devices, such as a PC, may not be able to continuously stream music	When using your Bluetooth phone as music player, disconnect it from other devices, such as PCs

Cause	Possible remedy
The volume of the music is too loud or too quiet. The volume of phone conversations is not comfortable (too low or too high)	
Different audio sources have different outputs	Reduce the hearing aid volume by pressing the (–) button. If the sound volume over the phone (or music) is too low after this reduction use the phone volume button to increase the volume. To change the volume of the phone or any other sound source streaming to uDirect 3, first adjust the overall level of the device and make minor adjustments on your uDirect 3
I have difficulties understanding the phone conversation in a noisy environment	
Your hearing aid microphones are set too high during the call	Reduce the environmental noise by pressing the volume down button for two seconds and then use the phone volume button to increase the phone volume

Cause	Possible remedy
The caller can hear me but I cannot hear the caller	
The call was transferred to your phone	Make sure the call is transferred back to uDirect 3 via the phone menu or long press of (A) button
uDirect 3 is no longer recognized by the phone or other Bluetooth device to which it was previously paired	
Two minutes after no device was found, uDirect 3 stops seeking devices to save power	Switch uDirect 3 off, wait until the red indicator switched off, then turn device on again
The devices are too far away from each other	Bring devices within one meter of each other and try again
uDirect 3 can be paired to a maximum of eight different Bluetooth devices. If the internal memory is full, new pairings may overwrite previous pairings	Repeat the pairing between uDirect 3 and the Bluetooth device which was over-written
Pairing was deleted	Repeat the pairing between uDirect 3 and the Bluetooth device which was deleted

Cause	Possible remedy
I am not streaming but the hearing aids keeps switching between the Bluetooth and the normal hearing aid program	
Your phone sends sounds to uDirect 3 due to incoming SMS or mails	Disable the system sounds for reminders, alarms, SMS, etc. via your phone menu
Phone key tones are switched on	Switch off phone key tones via your phone menu
I hear my phone's key-press tones in my hearing aids	
The key-press tones may be enabled on your phone	Deactivate key-press and confirmation tones on your phone. Please refer to the phone's user guide
Every time I switch on my phone, the music player starts	
For some phones this is normal behavior and not initiated by uDirect 3	Some phones can be configured not to start the music player automatically. Press and hold the access button to stop streaming. Stop the music player via phone menu after the phone is connected

Cause	Possible remedy
When I get into my car, the uDirect 3 connection to my phone shows unexpected behavior	
The car Bluetooth system connects to your phone	Consider disconnecting the car system from your phone
I accepted the incoming phone call but cannot hear the call in my hearing aids	
Some phones don't use uDirect 3 when accepting the call through the phone button	Always accept the call by pressing the uDirect 3 access (A) button

Cause	Possible remedy
The phone is ringing but the audio indicator is NOT solid blue and there is no ringing signal through the hearing aids	
Bluetooth is disabled in the phone	Ensure Bluetooth is enabled according to your phone's user guide
uDirect 3 and your phone are not connected to each other	1. Disable the Bluetooth functionality in your phone and turn it on again 2. Restart uDirect 3 by switching it off and then on. The uDirect 3 Bluetooth functionality will be automatically reactivated 3. Activate uDirect 3 by selecting it in the phone's "Active device" menu 4. Reduce the distance between uDirect 3 and the phone
uDirect 3 is not paired to the phone	Follow the pairing process in this guide
The distance is too great between uDirect 3 and the phone	Wear uDirect 3 around your neck. Keep distance to phone within <5m (15ft)

Cause	Possible remedy
The phone is ringing, the audio indicator is solid blue but there is no ringing signal through the hearing aids	
uDirect 3 is out of range of the hearing aids	Wear uDirect 3 correctly around your neck with both neckloop plugs firmly plugged in
Phone is set to vibrate	Enable ring tone on phone
If you are already on a phone call, then a second call cannot be recognized	Hang up the first phone call, then accept the second call
My phone rings but the name of the caller is not being spoken by uDirect 3	
The caller identification feature is not available or not enabled on your phone	Please refer to the phone's user guide to check the availability of this feature
Two phones are connected to your uDirect 3	If two phones are connected, caller identification is only available on the phone that was last paired to uDirect 3
Caller identification may be switched off	Please ask your hearing healthcare professional to enable caller identification on your uDirect 3

Cause	Possible remedy
During a phone call the other party reports difficulties hearing me	
uDirect 3 microphone openings may be covered	Make sure the microphone openings are not covered by part of your body, clothing or dirt and debris
uDirect 3 may not be worn correctly	Do not turn uDirect 3 sideways and make sure the neckloop plugs point towards your mouth as you speak
uDirect 3 may rub against clothing	Reduce movement while you speak, or consider using the external microphone (optional)
The surrounding noise may be too loud	Although uDirect 3 features noise reduction technologies, very noisy environments should be avoided. Suggest to your callers that they increase the phone volume, or hold the uDirect 3 closer to your mouth

Cause	Possible remedy
uDirect 3 is unresponsive or does not show any indicator when switching ON	
uDirect 3 may be completely discharged	Charge your uDirect 3 for at least 1 hour
It may indicate a software problem	Unplug any connectors from uDirect 3. Switch uDirect 3 off and on again. Consider resetting uDirect 3

The operating time of your uDirect 3 is reduced substantially

The battery typically needs several charging cycles to achieve full performance	Charge for at least 3 hours the first time you charge uDirect 3. Let a new uDirect 3 discharge completely three times
---	--

I have interruptions while my audio streaming is active

uDirect 3 is too far away from the hearing aids while streaming	Wear uDirect 3 clipped to your clothing with a maximum 40 cm (16") distance to the hearing aids
Bluetooth signal is out of range	Keep the Bluetooth device within the operating range

Cause	Possible remedy
I cannot hear the demo sound	
Hearing aids are not worn correctly	Make sure your hearing aids are inserted correctly
uDirect 3 is not worn correctly or switched off	Wear uDirect 3 clipped to your clothing with a maximum 40 cm (16") distance to the hearing aids. Make sure uDirect 3 is charged and switched on

For any problems not listed in this user guide, please contact your hearing healthcare professional.

Warnings

- 🔊 Please read the information on the following pages before using your Unitron uDirect 3.

Hazard warnings

- ⚠ Keep this device out of reach of pets and children under three years.
- ⚠ Only use hearing aids that have been specially programmed for you by your hearing healthcare professional.
- ⚠ Changes or modifications to any of the devices that were not explicitly approved by Unitron are not permitted.
- ⚠ Do not wear uDirect 3 around the neck while it is connected to any programming equipment.
- ⚠ Make sure to always remove both plugs of the neckloop and other cables when using the uDirect 3 as a handheld remote control.

- ⚠ Opening the uDirect 3 might damage it. If problems occur which cannot be resolved by following the remedy guidelines in the troubleshooting section of this user guide, consult your hearing healthcare professional.
- ⚠ This product has an embedded, non-replaceable battery. Do not attempt to open the product or remove the battery as this may cause injury and damage the product.
- ⚠ Dispose of electrical components in accordance with your local regulations.
- ⚠ Do not use the device in explosive areas (mines or industrial areas with danger of explosions, oxygen rich environments or areas where flammable anesthetics are handled) or where electronic equipment is prohibited.
- ⚠ External devices may only be connected if they have been tested in accordance with corresponding IECXXXXX standards. Only use accessories approved by Unitron.

- ⚠ Using your uDirect 3 cables in any way contradictory to their intended purpose (e.g., wearing the USB cable around the neck) can cause injury.
- ⚠ When operating machinery, ensure that no parts get caught in the machine.
- ⚠ The mini-USB port is to be used for the described purpose only.
- ⚠ Caution: electric shock. Do not insert plug alone into electrical outlets.

- ⚠ For safety reasons, only use chargers supplied by Unitron or stabilized chargers with a rating of 5VDC, min. 500 mA.

Information on product safety

- ⚠ Check your hearing aid compatibility with your hearing healthcare professional.
- ⚠ Protect the uDirect 3 from excessive moisture (bathing, swimming), heat (radiator, car dashboard), and direct skin contact when sweating (workout, fitness, sports).
- ⚠ Special medical or dental examination including radiation described below, may adversely affect the correct functioning of your remote. Remove and keep it outside the examination room/area before undergoing a medical or dental examination with X-ray (also CT scan) or any medical examinations with MRI/NMRI scans, generating magnetic fields.
- ⚠ Do not drop the device. Dropping onto a hard surface can damage your device.
- ⚠ Protect the uDirect 3 from excessive shock and vibration.

- ⚠ Do not use excessive force when connecting your uDirect 3 to the different cables.
- ⚠ Protect all openings (microphones, antenna, audio, FM and charger) from dirt and debris.
- ⚠ Never use a microwave or other heating devices to dry the uDirect 3.
- ⚠ Clean the uDirect 3 using a damp cloth. Never use household cleaning products (washing powder, soap, etc.) or alcohol to clean the uDirect 3.
- ⚠ The digitally-coded, inductive transmission technology used in the uDirect 3 is highly reliable and experiences virtually no interference from other devices. It should be noted, however, that when operating the hearing system near computer equipment, larger electronic installation or other strong electromagnetic fields, it may be necessary to be at least 60 cm (24") away from the interfering device to ensure proper operation. Keep the device at least 10 cm away from any kind of magnets.

68

- ⚠ Do not disconnect the neckloop while the uDirect 3 is transmitting signals to your hearing aids.
- ⚠ For safety reasons, recharge the uDirect 3 only with chargers supplied by Unitron or with USB-certified chargers >500 mA.
- ⚠ Do not connect an USB or audio cable exceeding 3 meters (9ft) in length to the uDirect 3.
- ⚠ When the uDirect 3 is not in use, turn it off and store it safely.

Other important information

- ⚠ High-powered electronic equipment, larger electronic installations and metallic structures may impair and significantly reduce the operating range.
- ⚠ If the hearing aids do not respond to the uDirect 3 because of an unusual field disturbance, move away from the disturbing field.

69

- ⚠ Your hearing aids and uDirect 3 are given a unique communication code during the fitting. This ensures that the device will not affect hearing aids worn by others.
- ⚠ When using an FM transmitter, be aware that radio signals might also be picked up and overheard by other receivers.

Information and explanation of symbols

With the CE symbol, Unitron confirms that this Unitron product – including accessories – meets the requirements of the Medical Devices Directive 93/42/EEC as well as the Radio Equipment Directive 2014/53/EU on radio and telecommunications transmitters.

This symbol indicates that it is important for the user to read and take into account the relevant information in this user guide.

This symbol indicates that it is important for the user to pay attention to the relevant warning notices in this user guide.

Important information for handling and product safety.

Australia and New Zealand EMC and Radiocommunications compliance label.

This symbol shall be accompanied by the name and the address of the authorised representative in the European Community.

This symbol shall be accompanied by the name and the address of the manufacturer (who are placing this device on the market).

Operating conditions:

This device is designed such that it functions without problems or restrictions if used as intended, unless otherwise noted in these user guides.

0° to 45° Celsius and relative humidity of <95% (non condensing).

Transport and storage conditions:

During transport or storage, the temperature should not exceed the limit values of -20°/60° Celsius and relative humidity of 90% for a long period of time. The air pressure between 500 and 1100 hPa is appropriate.

The symbol with the crossed-out garbage bin is to make you aware that this device may not be thrown away as normal household waste. Please dispose of old or unused devices, at waste disposal sites intended for electronic waste, or give your device to your hearing healthcare professional for disposal. Proper disposal protects the environment and health.

The Bluetooth® word mark and logos are registered trademarks owned by Bluetooth SIG, Inc. and any use of such marks by Unitron is under license. Other trademarks and trade names are those of their respective owners.

Apple and the Apple logo are trademarks of Apple Inc., registered in the U.S. and other countries. App Store is a service mark of Apple Inc.

Google Play is a trademark of Google Inc.

Compliance information

Declaration of Conformity

Hereby, Unitron Hearing declares that this product meets the requirements of the Medical Device Directive 93/42/EEC and complies with Radio Equipment Directive 2014/53/EU. The full text of the Declarations of Conformity can be obtained from the manufacturer.

Notice 1:

This device complies with Part 15 of the FCC Rules and with RSS-210 of Industry Canada. Operation is subject to the following two conditions:

- 1) this device may not cause harmful interference, and
- 2) this device must accept any interference received, including interference that may cause undesired operation.

Notice 2:

Changes or modifications made to this device not expressly approved by Unitron may void the FCC authorization to operate this device.

Notice 3:

This device has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules and ICES-003 of Industry Canada. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This device generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this device does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

Reorient or relocate the receiving antenna.

Increase the separation between the device and receiver.

Connect the device into an outlet on a circuit different from that to which the receiver is connected.

Consult the dealer or an experienced radio/TV technician for help.

Australia: Supplier Code Number N15398
New Zealand: Supplier Code Number Z1285

Unitron Hearing, a division of
National Hearing Services Inc.
20 Beasley Drive, P.O. Box 9017,
Kitchener, ON N2G 4X1 Canada

Sonova AG
Laubisrütistrasse 28 CH-8712
Stäfa, Switzerland

For a listing of Unitron group
companies, please visit
www.unitron.com

Appendix 3 – 029-5990-02b_uDirect3_Gd_ENG

Distributor

029-5990-02b_uDirect3_Gd_ENG © 2015, 2016 Unitron. All rights reserved.

unitron.com

Appendix 4 – 029-5992-02d_uTV3_Gd_ENG

uTV3 guide

Amir Panah

unitron. Hearing matters

Thank you

Thank you for choosing the Unitron uTV~ 3.

At Unitron, we care deeply about people with hearing loss. We work closely with hearing healthcare professionals to make advanced, purpose-driven solutions available to everyone. Because hearing matters.

The intended use of the uTV 3 is to connect your hearing aids to a TV or other audio source. It sends audio signals wirelessly to the hearing aids through a streamer within a distance of 30 m (100 ft). It also serves as a charging station for the uDirect~ 3 and uStream.

Your uTV 3

Hearing healthcare professional: _____

Telephone: _____

Serial number of uTV 3: _____

Warranty: _____

Date of purchase: _____

Your streamer:

☐ uStream

☐ uDirect 3

Package contents

The complete Unitron uTV 3 consists of the following elements:

uTV 3 plus the following cables:


- Power supply with mini-USB plug and country adapter
- “Optical”
- “Coaxial”
- “Analog RCA/Cinch”
- “Analog 3.5mm”

Table of contents

uTV 3 description	2
Getting started	4
Step 1: Setting up the power supply	4
Step 2: Connecting to the TV	6
Step 3: Charging your streamer	12
Step 4: Pairing uTV 3 with your streamer.....	15
Daily usage of the uTV 3	17
Step 1: Wearing your streamer	17
Step 2: Listening to the TV	18
Step 3: Adjusting the volume of the TV signal.....	20
Step 4: Switching off uTV 3 and streamer.....	22
Understanding the indicator lights	23
Testing the audio transmission	24
Troubleshooting guide	26
Important information.....	29
Information and explanation of symbols.....	33
Compliance information.....	35


uTV 3 description

The Unitron uTV 3 connects your hearing aids wirelessly to your TV or other audio devices through a streamer. The system consists of three main elements: your hearing aids, a streamer and the uTV 3.

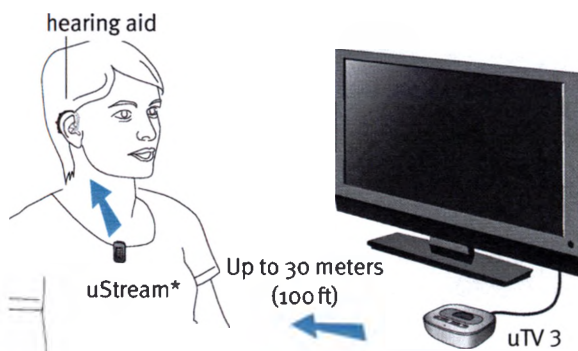*uDirect 3 is also compatible.

The uTV 3 is connected to your TV (or any other audio source) and wirelessly sends audio signals from your TV to a streamer within a range of up to 30 meters (100 ft). The uTV 3 also serves as a charging cradle for your streamer when it is not in use.

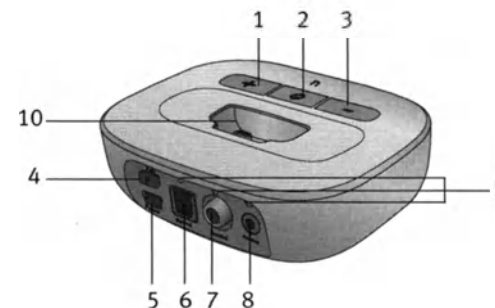- 1 Volume up (+) button
- 2 Power on/off button and main indicator light
- 3 Volume down (-) button
- 4 Pairing (P) button
- 5 Power supply
- 6 Digital audio (optical)
- 7 Digital audio (coaxial)
- 8 Analog audio (3.5 mm)
- 9 Audio indicator lights
- 10 Exchangable charging slot

Getting started

Step 1: Setting up the power supply

To attach the country adapter:

1. Select the adapter that corresponds to your country. (Included country adapters may vary by country.)

2. Insert the round corner of the adapter to the round corner of the universal power supply until it is completely inserted.
3. Click the tail of the adapter in to lock it into place. Make sure the adapter is securely locked.

To remove the country adapter:


1. Pull back the tab on the universal power supply and hold.
2. Gently pull the adapter upwards to remove it from the receptacle.

Connecting the power supply

Place the uTV 3 next to your TV.

Connect the mini-USB plug of the power supply to the power input (Power) of the uTV 3. Plug the power supply into the wall socket.

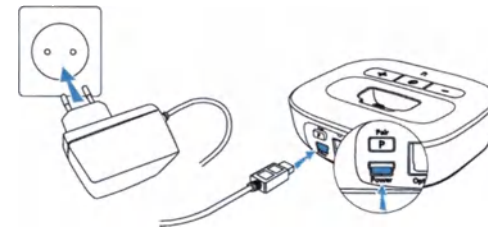Step 2: Connecting to the TV

Supported audio formats:

The uTV 3 can be connected to any audio source such as a TV, PC or hi-fi system. It supports three different audio input formats.


For best results, we recommend you use the optical or coaxial connection if available on your audio device.

 **Note:** please ensure to always connect uTV 3 to the audio output of your device.

Signal priorities:

Normally it is enough to connect one audio cable. If more than one cable is connected and a signal is present, uTV 3 will select the best possible input in the following order: optical, coaxial, analog.

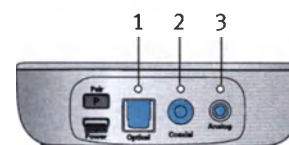 **Note:** even though the optical and coaxial cables are only one wire, the signal is stereo.


uTV 3 shows the selected input by an orange audio indicator light:

Input selected	Slow blink orange	
Input selected and audio available	Solid orange	

Option 1: connecting via optical cable

If an optical output is available:

- Unpack the cable labeled “Optical”
- Remove the protection caps from the plugs
- Connect one end of the optical cable into the “Optical” input of uTV 3, and the other end into the optical output of your TV

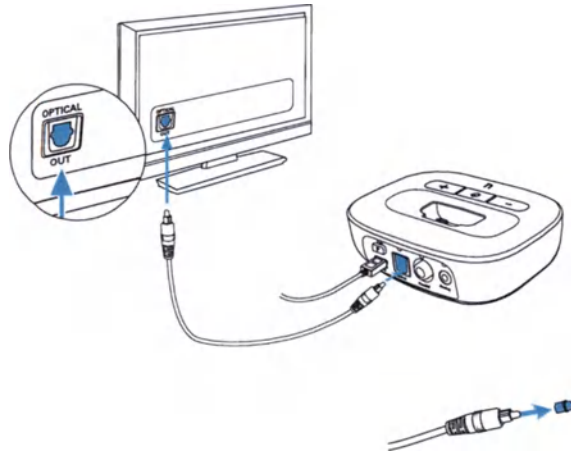
8

Option 2: connecting via coaxial cable

If a coaxial output is available:

- Unpack the cable labeled “Coaxial” (orange plugs)
- Connect one end of the coaxial cable into the “Coaxial” input of uTV 3, and the other end into the coaxial (also known as “digital out” or “S/P DIF”) output of your TV

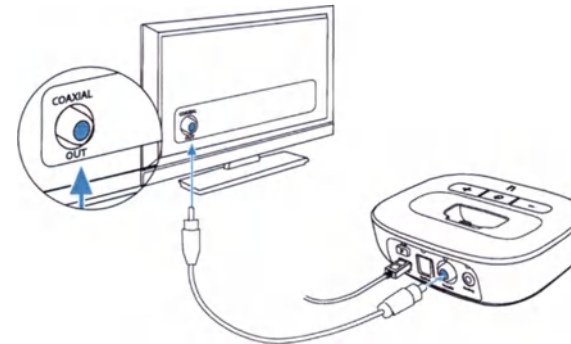
9

Option 3: connecting via RCA/Cinch cable

If an RCA output (red / white) is available:


- Unpack the cable labeled “Analog RCA/ Cinch”
- Connect the 3.5 mm plug into the “Analog” input of the uTV 3. Connect the red and white plugs to the matching red and white output sockets on your TV

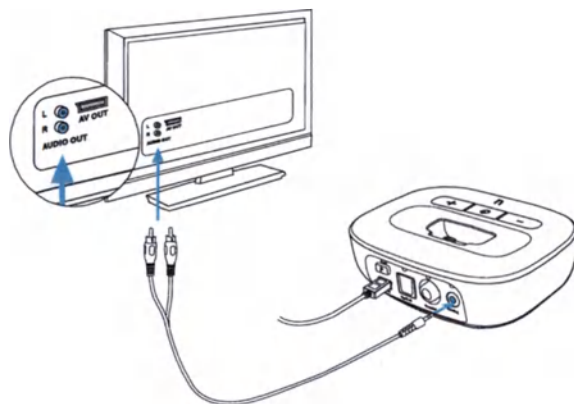10

Option 4: connecting to the headset socket

- Unpack the cable labeled “Analog 3.5 mm”.
- Connect the stereo plug into the “Analog” input of the uTV 3. Plug the other end into the headset output socket of your TV

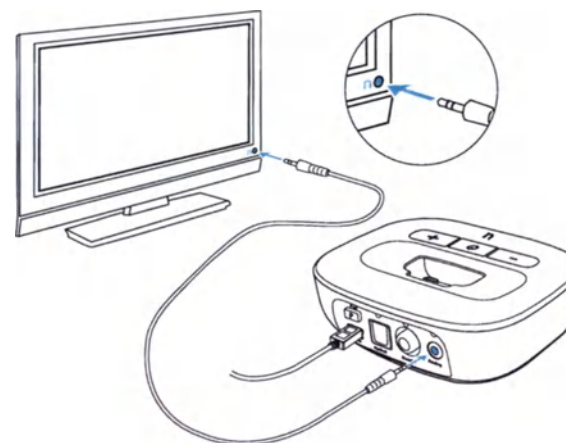11

-  When connecting the uTV 3 to the headset socket of your TV, you are able to adjust the volume via your TV remote control.
-  Some TVs will switch off their loudspeakers when plugging the jack into the headset socket – the TV will not be audible for other people. In this case, use the Cinch/RCA output socket as described in Option 3. For further information on volume adjustment, refer to the section on adjusting the volume. The uTV 3 can also be connected to any other audio source such as an MP3 player, PC or home sound system using the stereo plug cable.

Step 3: Charging your streamer


To charge the streamer, gently plug it into the charging slot of the uTV 3 as shown below. Ensure the integrated plug is inserted into the streamer charging socket correctly.

Ensure that you are using the correct charging slot.


Changing the charging slot for your streamer

The uTV 3 comes with two charging slots. You may need to change the slot using a Phillips head screwdriver.

1. Unplug the USB power cable.
2. Unscrew the screw on the bottom of the uTV 3. Lift the charging slot to remove.

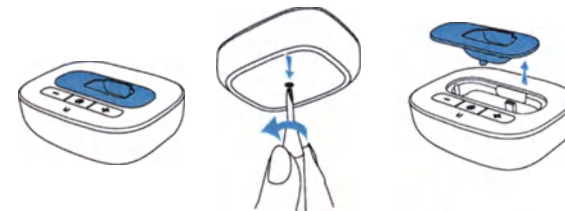3. **Insert the other slot into the opening.** Replace the screw at the bottom and gently tighten it.

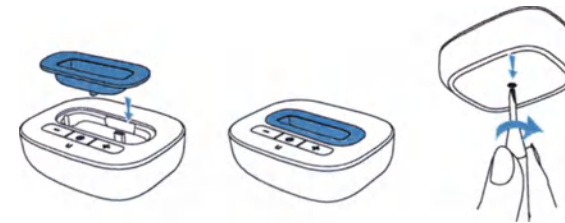-  When using your streamer for the first time, charge it for at least three hours even if the battery indicator turns green earlier. After the initial charge, you may leave the streamer in the charging slot for and unlimited time.

uStream

uDirect 3

Charging is complete when the streamer's indicator light turns solid green.

The uTV 3 can either be on or off during charging. Its status indicator light does not indicate the charging process.

The typical time to fully charge the streamer is 90 minutes.

We recommend that the streamer is switched off while charging.

-  For safety reasons, recharge the battery only with chargers supplied by Unitron. If other chargers are used they must be USB certified (5.0 VDC / 500mA).

Step 4: Pairing uTV 3 with your streamer


-  If you have received the uTV 3 and the Unitron streamer together in a complete set, they are paired already. If they were purchased separately, they must first be paired with each other.

First, put your streamer into pairing mode. (Check your streamer's user guide for further details.)

Next, put the uTV 3 into pairing mode.

To manually start pairing mode on the uTV 3:

1. Switch uTV 3 on.
2. Press the pairing button on the back, and the main indicator light will start rapidly blinking blue.
3. Make sure uTV 3 and the streamer are within 1 m (3 ft) range.

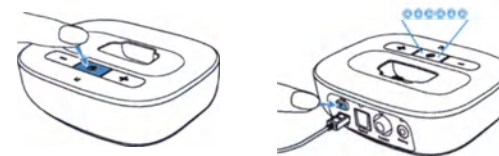- 1 uTV 3 will remain in pairing mode for two minutes or until the pairing has been completed.

Completion of pairing


Wait until both indicator lights turn solid blue. This can take up to two minutes but usually occurs within a few seconds.

- 1 The uTV 3 basestation can be paired to a maximum of four streaming devices. A fifth pairing will replace the first paired device.
- 1 The two devices only need to be paired once. They remain paired after switching off and will automatically connect after being switched on again.

Daily usage of the uTV 3

Step 1: Wearing your streamer

Wear the streamer as shown and switch it on.

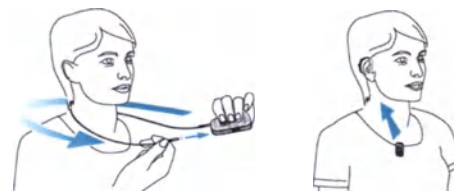- 1 For a more detailed description of the streamers, refer to the streamer's user guide.
- 1 After switching your streamer on, it will automatically recognize your Unitron wireless hearing aids.
- 1 In cases of interference with other hearing aids, contact your hearing healthcare professional.
- 1 The uTV 3 can only transmit the audio signal to one streamer at a time. Streaming to multiple receiving streamers is not possible.

Step 2: Listening to the TV

You will hear the sound of the TV directly in your hearing aids.

 If you cannot hear the sound of your TV through uTV 3, repeat the setup in steps 1-4.

 If no audio is available, uTV 3 automatically stops the Bluetooth audio transmission to the streamer after five minutes. Its main indicator light will blink green.

The uTV 3 supports Bluetooth operation up to a range of 30 m (100ft). Line-of-sight between your streamer and the uTV 3 is not required. However, interference in the environment by walls or furniture may reduce the operating range.

If the streamer moves out of range of the uTV 3, the wireless audio transmission will be terminated. If the streamer is then moved back into range of the uTV 3, it will automatically resume streaming to your hearing aids.

 It may take up to ten seconds until the streaming resumes. Both streaming light indicators, should be solid blue during streaming.


 If the uTV 3 status light:

- blinks green, make sure the audio input to the uTV 3 is properly connected.
- is green, make sure the streamer is switched on, both ends of the neckloop are plugged in (for uDirect 3 or uDirect 2) and the devices are paired.
- blinks blue, you are in pause mode. Press the streamer access (A) button to resume streaming.

Step 3: Adjusting the volume of the TV signal

First, set the volume of the TV speakers at the level to suit other viewers.

Next, adjust the volume via the uTV 3 volume control until the sound through the hearing aids is comfortable.

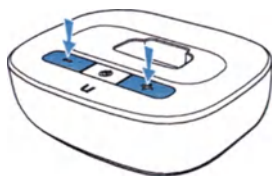To fine-tune the overall volume of the hearing aids and uTV 3, you may use the streamer's volume (+/-) buttons.

20

If your uTV 3 is connected to the headphone socket of your TV, you may also use your TV remote control to adjust volume via the TV directly. Refer to the TV user guide for more details on setting the TV volume.


-  If your TV is connected to a satellite receiver or a set-top box, make sure the volumes of these additional devices are set to a 75% level for best audio quality results.
-  Refer to your streamer's user guide for pairing instructions of the streamer with a mobile phone or other entertainment and communication devices such as radios or personal computers.
-  Refer to your streamer's user guide for information on using the streamer to pause the audio transmission and other functions.

21

Step 4: Switching off uTV 3 and streamer

To switch off the devices:

To end the transmission of the TV audio signal, switch off the uTV 3 by pressing the on/off button as shown .

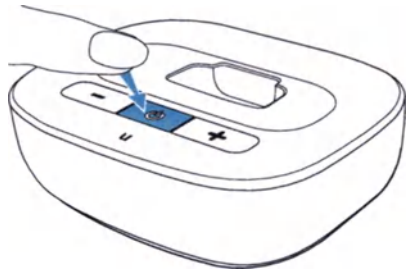Understanding the indicator lights

Switched on	Solid green	
On but no audio input	Short blink green	
Pairing mode	Rapid blue blinking	
Streaming mode	Solid blue	
Pause mode	Short blink blue	
Test mode	Multi-color blink	

For a more detailed explanation of the streamer indicator lights, please refer to the streamer's user guide.

Testing the audio transmission

The uTV 3 is able to generate a test sound to check the system.

To start test mode:

- Switch on and wear your streamer correctly.
- Switch off the uTV 3 first, then press the uTV 3 volume (—) and power (🔌) button simultaneously until the status light starts to blink in different colors.
- The streamer and the uTV 3 should connect now (blue/green or blue/orange blinks) and you should hear a sound test melody through your hearing aids.
- The test sound continues for two minutes or until any of the uTV 3 buttons are pressed.
- Use the streamer volume (+/–) buttons to change the volume of the test sound.
- If no test sound can be heard, redo the pairing procedure.

To test the cable connections:

If you were able to hear the test melody from uTV 3 but you cannot hear the sound from the TV, there is a problem with the cable connection.

- Ensure the audio cable is completely plugged in
- Ensure that you connected the cable to an output and not an input
- Try a different cable
- Try a different output socket of the TV

Color codes during test mode:

Red: no audio input on uTV 3

Green: audio input detected on uTV 3

Orange: searching for a streamer

Blue: streaming test sound to streamer

 If no streamer is found and no audio input is available the uTV 3 will stop the test mode and its light indicator will blink green.

Troubleshooting guide

Most probable cause	Remedy
None of the indicator lights turned solid blue	
Devices are turned off	Turn on devices
Devices are not paired	Repeat pairing
I see a rapidly blinking blue indicator light on the uTV 3	
uTV 3 is in Bluetooth pairing mode	Go to the pairing section of your streamer's user guide
Only one of the blue indicator lights turned solid blue	
One of the devices is connected to another Bluetooth device	Make sure to switch off all other Bluetooth devices except the streamer and uTV 3
Both indicator lights are solid blue but I cannot hear the TV	
The streamer is out of range of the hearing aids	Wear your streamer correctly. Turn all devices off and on again
The streamer is not correctly programmed for your hearing aids	Ask your hearing healthcare professional to regroup your streamer and hearing aids

Most probable cause	Remedy
I cannot hear the sound from the TV via the streamer	
Devices are not on	Turn on devices
Cables not plugged in correctly	Repeat steps 2 and 3
No audio input to uTV 3, the light indicator blinks green	Make sure TV is not muted Select a 75% volume level on the TV Make sure the audio cable uses an "out" port
uTV 3 is in pause mode	Indicated by blue blinking status lights on both devices. A short press of the streamer access (A) button will resume streaming
Distance between streamer and uTV 3 is too far	Bring devices closer to each other, 30 m (100 ft) is the maximum in line of sight To avoid streaming interruptions, make sure you wear your streamer correctly and keep the uTV 3 within the operating range
Some TVs may not output stereo sound when in surround sound mode	Try switching TV or devices to stereo output

Most probable cause	Remedy
---------------------	--------

Volume during streaming is not comfortable
(too low or too high)

TV audio sources have different output volumes	If uTV 3 is connected to the optical, coaxial or RCA/ Cinch socket of your TV, adjust the volume by using the uTV 3 volume buttons If the uTV 3 is connected to the headset output of the TV, you may adjust the volume using the TV volume control Adjust volume of your satellite receiver or set-top box Adjust the volume using the streamer volume buttons
--	---

After connecting the uTV 3, the TV loudspeaker is switched off

The headset output connection of the TV is selected	Choose another audio output on your TV such as optical, coaxial or RCA/ Cinch
---	---

Important information

- ⚠ Please carefully read the corresponding section in the uDirect 3 or uStream user guide which is delivered along with each accessory.

Hazard warnings

- ⚠ Keep this device out of reach of pets and children under three years.
- ⚠ Only use hearing aids that have been specially programmed for you by your hearing healthcare professional.
- ⚠ Changes or modifications to any of the devices that were not explicitly approved by Unitron are not permitted.
- ⚠ Opening the uTV 3 might damage it. If problems occur which cannot be resolved by following the remedy guidelines in the troubleshooting section of this user guide, consult your hearing healthcare professional.

- ⚠ Dispose of electrical components in accordance with your local regulations.
- ⚠ Do not use the device in explosive areas (mines or industrial areas with danger of explosions, oxygen rich environments or areas where flammable anesthetics are handled) or where electronic equipment is prohibited.
- ⚠ External devices may only be connected if they have been tested in accordance with corresponding IECXXXXX standards. Only use accessories approved by Unitron.
- ⚠ Using your uTV 3 cables in any way contradictory to their intended purpose (e.g., wearing the USB cable around the neck) can cause injury.
- ⚠ When operating machinery, ensure that no parts get caught in the machine.
- ⚠ The mini-USB port is to be used for the described purpose only.

- ⚠ Caution: electric shock. Do not insert plug alone into electrical outlets.

- ⚠ For safety reasons, only use chargers supplied by Unitron or stabilized chargers with a rating of 5VDC, min. 500 mA.

Information on product safety

- ⚠ When uTV 3 is not in use, turn off the device.
- ⚠ Always plug in the uDirect 3 neckloop and wear it around your neck when transmitting signals to your hearing aids.
- ⚠ Protect the device connectors, plugs, charging slot and power supply from dirt and debris.
- ⚠ Protect all devices from excessive moisture (bathing or swimming areas) and heat sources (radiator, TV top). Protect devices from excessive shock and vibration.

- ⚠ Clean devices using a damp cloth. Never use household cleaning products (washing powder, soap, etc.) or alcohol to clean them. Never use a microwave or other heating devices to dry any of the devices.
- ⚠ The digitally-coded, inductive transmission technology used in these devices is extremely reliable and experiences virtually no interference from other devices. It should be noted, however, that when operating the system near computer equipment, larger electronic installation or other strong electromagnetic fields, it may be necessary to be at least 60 cm (24") away from the interfering device to ensure proper operation.
- ⚠ X-ray radiation, CT or MRI scans may destroy or adversely affect the correct functioning of the devices.
- ⚠ High-powered electronic equipment, larger electronic installations and metallic structures may impair and significantly reduce the operating range.

Information and explanation of symbols

With the CE symbol, Unitron confirms that this Unitron product – including accessories – meets the requirements of the Medical Devices Directive 93/42/EEC as well as the Radio Equipment Directive 2014/53/EU on radio and telecommunications transmitters.

The numbers after the CE symbol correspond to the code of certified institutions that were consulted under the above-mentioned directives.

This symbol indicates that it is important for the user to read and take into account the relevant information in this user guide.

This symbol indicates that it is important for the user to pay attention to the relevant warning notices in this user guide.

Important information for handling and product safety.

Australia and New Zealand EMC and Radiocommunications compliance label.

This symbol indicates that the products described in these user instructions adhere to the requirements for an application part of Type B of EN 60601-1. The surface of the hearing aid is specified as applicated part of Type B.

This symbol shall be accompanied by the name and the address of the authorised representative in the European Community.

This symbol shall be accompanied by the name and the address of the manufacturer (who are placing this device on the market).

Operating conditions:

This device is designed such that it functions without problems or restrictions if used as intended, unless otherwise noted in these user guides.

Transport conditions:

Temperature: -20° to $+60^{\circ}$ Celsius (-4° to $+140^{\circ}$ Fahrenheit).

Humidity: Up to 90% (non condensing).

The symbol with the crossed-out garbage bin is to make you aware that this device may not be thrown away as normal household waste. Please dispose of old or unused devices, at waste disposal sites intended for electronic waste, or give your device to your hearing healthcare professional for disposal. Proper disposal protects the environment and health.

The Bluetooth® word mark and logos are registered trademarks owned by Bluetooth SIG, Inc. and any use of such marks by Unitron is under license. Other trademarks and trade names are those of their respective owners.

Compliance information

Declaration of Conformity

Hereby, Unitron Hearing declares that this product meets the requirements of the Medical Device Directive 93/42/EEC and complies with Radio Equipment Directive 2014/53/EU. The full text of the Declarations of Conformity can be obtained from the manufacturer.

Notice 1:

This device complies with Part 15 of the FCC Rules and with RSS-210 of Industry Canada. Operation is subject to the following two conditions:

- 1) this device may not cause harmful interference, and
- 2) this device must accept any interference received, including interference that may cause undesired operation.

Notice 2:

Changes or modifications made to this device not expressly approved by Unitron may void the FCC authorization to operate this device.

Notice 3:

This device has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules and ICES-003 of Industry Canada.

These limits are designed to provide reasonable protection against harmful interference in a residential installation. This device generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions,

Appendix 4 – 029-5992-02d_uTV3_Gd_ENG

may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this device does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.

- Increase the separation between the device and receiver.

- Connect the device into an outlet on a circuit different from that to which the receiver is connected.

- Consult the dealer or an experienced radio/TV technician for help.

Notice 4:

FCC/IC RF Radiation Exposure Statement

This equipment complies with FCC/Industry Canada RF radiation exposure limits set forth for an uncontrolled environment.

This equipment should be installed and operated with a minimum distance of 20 cm between the radiator and your body.

This transmitter must not be co-located or operating in conjunction with any other antenna or transmitter.

Australia: Supplier Code Number N15398

New Zealand: Supplier Code Number Z1285

Unitron Hearing, a division of
National Hearing Services Inc.
20 Beasley Drive, P.O. Box 9017,
Kitchener, ON N2G 4X1 Canada

Sonova AG
Laubisrütistrasse 28
CH-8712 Stäfa, Switzerland

For a listing of Unitron group
companies, please visit
www.unitron.com

Appendix 4 – 029-5992-02d_uTV3_Gd_ENG

Distributor

D/2017-05 029-5992-02 © 2015-2017 Unitron. All rights reserved.

unitron.com

Appendix 5 – 029-5993-02c_uStream_Gd_ENG

uStream
guide

Amir Panah

unitron. Hearing matters

Thank you

Thank you for choosing the Unitron uStream.
At Unitron, we care deeply about people with hearing loss. We work closely with hearing healthcare professionals to make advanced, purpose-driven solutions available to everyone. Because hearing matters.

The intended use of uStream is to receive inputs from different Bluetooth® devices and transmit them wirelessly to your hearing aids.

Virtually any Bluetooth enabled audio device such as mobile phones or MP3-players are supported. Enjoy speech and music from other multimedia devices, such as computers, TV or music players. uStream can also be used as a remote control to adapt your hearing aids to different listening situations.

 Please note that your hearing healthcare professional must configure your uStream before it can be used as a remote control for your hearing aids.

Your uStream

Hearing healthcare professional: _____

Telephone: _____

Serial number of uStream: _____

Warranty: _____


Date of purchase: _____

The uControl™ app can be downloaded at the App Store and Google play. Search for Unitron and select uControl.

Quick reference

Charging uStream

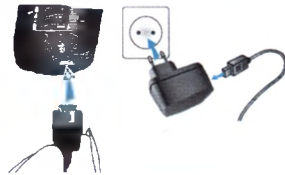
Wearing uStream

Bluetooth pairing

On/Off

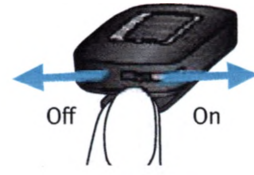Indicator light

- Solid red - charging
- Solid green - fully charged
- Red blink - low battery
- No light - battery drained
- Solid pink - demo sound
- Solid blue - streaming
- Blue short blink - 1 device connected
- Blue double blink - 2 devices connected
- Blue rapid blink - pairing

Table of contents

Your uStream at a glance	2
Setting up the power supply	4
Charging the uStream	5
Turning the uStream on and off	6
Summary of indicator lights	7
Wearing the uStream correctly.....	9
Operating distance.....	9
Button lock function.....	10
Testing sound via the uStream.....	11
Access button functionality	12
Spoken notifications	12
uControl app.....	13
Using uStream as a remote control for your hearing aid	14
Using uStream with a mobile phone	16
Advanced phone options	23
Using uStream with an audio device.....	27
Using uStream with uTV 3	30
Using uStream with uMic	33
Managing multiple audio inputs to the uStream.....	36
Deleting the pairing of Bluetooth devices	38
Troubleshooting guide	39
Warnings	47
Information and explanation of symbols.....	53
Compliance information.....	55
Patient feedback	57
Additional notes	58

Your uStream at a glance

Buttons


- 1 Access (A)
- 2 Volume up/down (+/-)
- 3 On/Off switch (VO)

Inputs

- 4 Microphone openings
- 5 Charging input (mini-USB)


Others

- 6 Indicator light
- 7 Eyelet for lanyard
- 8 Retaining clip

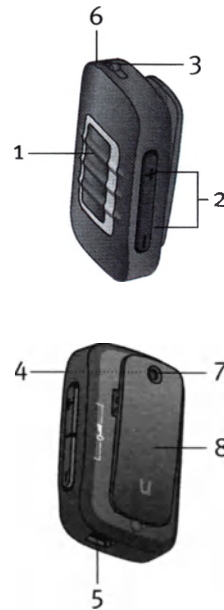Accessories

- 9, 10 Universal power supply with mini-USB plug
- 11-14 US, EU, UK, AUS adapter

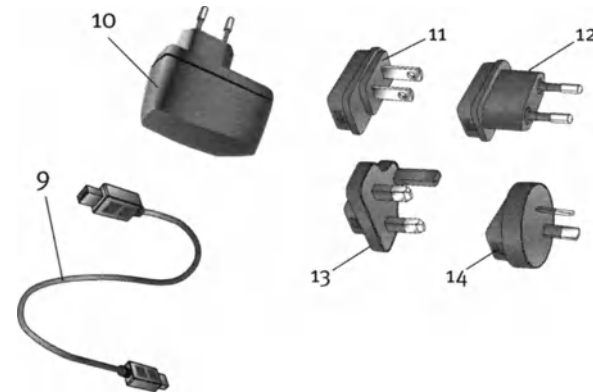
Included country adaptors may vary by country.

Setting up the power supply

Select the adapter that corresponds to your country.


1. Insert the round corner of the adapter to the round corner of the universal power supply until it is completely inserted.
2. Click the tail of the adapter in to lock it into place. Make sure the adapter is securely locked.

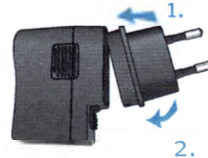To remove the power adapter:

1. Pull back the tab on the universal power supply and hold.
2. Gently pull the adapter upwards to remove it.

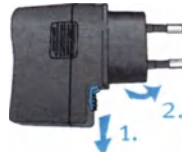Charging the uStream

Your uStream has a non-removable rechargeable battery.

To charge it:

1. Plug the larger end of the charging cable into the universal power supply. Plug the smaller end plug into the mini-USB port on uStream.

2. Plug the power supply into a power outlet.

During charging, the indicator light is red and when the battery is fully charged, the indicator light turns green.

It is safe to leave the uStream connected to the charger overnight.

When charging your uStream for the first time, let it charge for at least three hours.

uTV™ 3 can also be used as a charging station for uStream.

Turning the uStream on and off

Slide the power switch (I/O) as shown to turn the uStream on or off.

Switch the uStream off when not in use or if in storage or during transportation.

For operational use switch the uStream on.

Summary of indicator lights

The indicator light provides helpful information as follows:

Indicator light

Charging	Solid red	
Charged and full battery	Solid green	
Switching on	Green 2 seconds	
Switching off	Solid red two seconds	
Battery above 10%, normal operation	Slow blink green	
<10% battery, charging required (at least two hours)	Short blink red**	
Battery dead, charging required (at least two hours)	No light after switching on	
Enable/disable lock state	Short blink red (lock) or short blink green (unlocked)	
Bluetooth phone or music streaming	Solid blue	

One Bluetooth device connected	Short blink blue	
Two Bluetooth devices connected	Double-blink blue	
Bluetooth pairing	Rapid blue blinking	
Demo sound playing	Solid pink	

Battery state

The color of the power indicator informs you about the uStream battery state.

The green indicator light indicates a battery level >10%.

The red indicator light indicates a battery level <10%. Less than one and a half hours streaming remain and the battery should be charged for at least two hours.

Wearing the uStream correctly

Before any phone, music or TV streaming can take place, insert your hearing aids into your ears and turn them on.

1. Switch uStream on.
2. Clip it securely onto your clothing. For best results, keep the uStream 20 cm (8") from your mouth. Do not exceed a distance of 40 cm (16").

 To prevent scratching sounds, do not cover uStream with clothing during phone calls.

Operating distance

- Optimum distance 20 cm/8"
- Maximum distance 40 cm/16"

Button lock function

The buttons on your uStream can be locked to avoid inadvertent presses. Press and hold the volume up and volume down (+ -) buttons simultaneously to lock or unlock the buttons. The indicator light will flash red three times to indicate that the buttons have been locked.

Press and hold the volume up and volume down (+ -) buttons simultaneously again to unlock the buttons. The indicator light will flash green three times to indicate that the buttons have been unlocked.

When locked, the indicator light will triple blink and you may hear a triple beep through the hearing aids each time a button is pressed.

You will still be able to answer an incoming telephone call using the access (A) button while it is locked.

 In a locked state, the indicator light will stop blinking.

 The key symbol on the back of uStream reminds you of the lock function.

Testing sound via the uStream

1. Wear the uStream correctly.
2. Press and hold the access (A) button while switching on the device. Release the access button. The indicator light will turn pink.

The hearing aids will switch to the streaming mode and you will hear the uStream sound demo.

Use the volume (+ -) buttons to adjust the volume.

The sound demo will stop after three minutes or when the access button is pressed. Adjusting the volume during the sound demo will also set the volume level of the spoken uStream messages.

If you cannot hear the sound demo, please see the troubleshooting guide.

Access button functionality

The access button has multiple functions:

- **Phone:** Manages telephone functions such as accepting and ending phone calls. See the section “Using uStream with a smartphone or cell phone” for more information.
- **Streaming:** Pauses and resumes the streaming from an audio device.
- **Remote control:** can be used to select the hearing aid programs.

Spoken notifications

Spoken notifications inform you about the different operation modes and status of the uStream using text-to-speech, for example low battery.

Your hearing healthcare professional can configure the language setting as well as which notifications are spoken.

12

 Caller identification is independent of the spoken notifications level.

 Spoken notifications must be enabled by your hearing healthcare professional.

uControl app

The uControl is a smartphone app which offers quick and easy volume and program adjustments.

It also provides:

- individual left and right volume control
- additional program controls
- easy Bluetooth pairing

13

Using uStream as a remote control for your hearing aid

 Before the remote control functionality can be used, your hearing healthcare professional must enable the remote control.

If uStream is used as a handheld remote control, the maximum distance to the hearing aids is 55 cm/22".
Hold the uStream as shown.

Volume change

Press the (+/−) volume buttons on the uStream to adjust the volume of your hearing aids:

- Press (+) to increase volume
- Press (−) to decrease volume

 Adjusting the volume on the uStream will have the same effect as adjusting the volume of the hearing aids.

Program change

Use the access (A) button to change hearing aid listening programs.

Changing hearing aid programs is only possible if no streaming is in progress.

More remote control information

- Your hearing aids may beep to confirm that your selection was applied.
- Both hearing aids are controlled at the same time.
- The available hearing aid programs and volume levels are set by your hearing healthcare professional.

Using uStream with a mobile phone

If you have a mobile phone equipped with Bluetooth functionality (e.g. iPhone, Android phone) uStream may be used to stream phone conversations wirelessly to your hearing aids allowing you to hear the callers voice directly in your hearing aids. uStream picks up your voice through its own microphone and streams it to your phone.

During phone calls, the access (A) button provides a range of features that help you manage your calls.

16

Pairing your mobile phone

1. Enable Bluetooth on your phone settings menu.
2. Starting with the uStream in the off position, press and hold the (+) button and turn the uStream on. The blue audio indicator will flash.
3. Enable Bluetooth pairing mode on your phone so that a new device can be added. See manufacturer instructions for details. If prompted, enter "0000" for pairing code.
4. If prompted, accept request to grant uStream access to your phone book. This is used for voice calling and caller ID features.
5. After successful pairing, the blue audio indicator on the uStream will stop blinking.

Once pairing is done, your phone and uStream should pair automatically.

 Please contact your hearing healthcare professional if you can not complete the pairing procedure with the steps described.

17

-  To receive phone calls your phone must support either the Bluetooth “Head Set” (HSP) or “Hands Free” (HFP) profile. When pairing to your phone, it may ask you which Bluetooth service you want to enable. Select “Headset” or, if available, “Stereo Headset.”
-  To access caller ID features, your phone must support the Phone Book Access Profile (PBAP)
-  When pairing to a cellular phone, we recommend disabling the phone key tones and SMS alerts.
-  To listen to stereo music from your phone, it must support the A2DP/“Stereo Headset” profile.

Preparing for phone calls

During phone calls with uStream, ensure your phone is within 5-10 m (15-30 ft) to ensure a reliable connection and ensure the uStream is clipped securely to your clothing within 40 cm of your hearing aids.

Making phone calls

Use the phone keypad on your phone to enter the number and press the phone dial button to initiate the call. Depending on the phone, you may hear the dial tone.

Speak into the microphone on the uStream – not the phone. In noisy environments, you may find that moving the uStream closer to your mouth will make it easier for the caller to hear your voice.

Make minor volume adjustments with the uStream volume buttons (+, -).

Incoming calls

The uStream must be switched on, paired to the phone and worn correctly.

When a call is received on your phone, the uStream audio indicator will show a solid blue light, and you will hear the ringtone signal through your hearing aids.

Accepting incoming calls

To accept an incoming call, briefly press the access (A) button on your uStream once you hear the ringtone through your hearing aids. It will take about five seconds before you hear the ringtone through your hearing aids.

Speak into the uStream, not into the phone, when the call is in progress.

-  For your safety, the hearing aid microphones will remain on at reduced levels during phone calls and audio streaming.
-  Only one phone call at a time is possible. A second ringing phone is ignored while another call is already connected.

Ending a call

To end a phone call via your uStream, briefly press the access (A) button. The call will terminate and your hearing aids will return to the previously-active program or audio source (such as the audio jack input).

You can also end a call by operating your phone instead of your uStream.

Different ringtones of the uStream

The phone ringtone depends on the uStream configuration and your phone capabilities.

The ringtone is assigned automatically during the initial phone pairing process and cannot be configured.

When pairing more than one phone, the second and third paired phone will have a different

ringtone. A total of three different ringtones are automatically assigned.

The uStream will also announce the caller name as shown on the phone display (caller identification using text-to-speech technology). The name is taken from your phone's directory.

 Please note that this feature must be supported by your phone. Please consult the phone user guide or the phone supplier for more information.

 The default language setting of a new uStream is US English. Your hearing healthcare professional may configure the appropriate language setting for proper pronunciation of the names. Your hearing healthcare professional may switch the caller identification off upon your request.

Rejecting a call

If you prefer not to answer an incoming call, press and hold the access (A) button on your uStream for two seconds until the ringtone stops. Using your phone to reject the call will have the same effect.

Advanced phone options

Your hearing healthcare professional may configure the following additional features on uStream. They will be assigned to a long press (two seconds) of the access (A):

 Your phone must support these features. Please consult the phone user manual for details.

Redial last number

A long press of the access (A) button will redial the last dialed number.

Voice dial

A long press of the access (A) button will start the voice dial which allows you to speak the person's name instead of typing the number.

 When using two phones at the same time, only the last paired phone can be used for voice dial.

Dial preferred number

A long press of the access (A) button will dial a pre-defined number.

Hold call

While in a phone conversation, a long press of the access (A) button will put the call on hold. The caller stays on the line but there is no audio transfer and the hearing aids switch to the previous hearing program. Use this mode to have a side conversation with a person next to you without transmitting the conversation to the person on the phone.

Transfer call

While in a phone conversation, a long press of the access (A) button will transfer the call to the phone. This allows you to hand your phone to another person. Your hearing aids will switch back to the previous hearing program.

After holding or transferring a call, another long press of the access (A) button transfers the call back to the uStream and the call is resumed through your hearing aids.

 As long as the call is on hold or transferred, you may use the access (A) button to change your hearing aid programs again.

Caller identification, spoken caller names

The uStream may be able to tell you the name of the caller.

If your phone supports Phone Book Access (PBA) or Phone Book Access Profile (PBAP), the uStream will try to access information about the caller every time the phone rings.

As the phone rings, the displayed information is transferred to the uStream and spoken through your hearing aids as follows:

- The stored name in your phone's phonebook. In the example, you will hear "Mike Office"
- The phone number, if the name is not known

To make sure the pronunciation is correct your hearing healthcare professional can configure the language settings of your uStream.

The phonebook is not downloaded to the uStream. This way, uStream can also support the use of different phones.

Please consult your phone's user guide or the phone dealer to check whether this feature is supported by your phone.

 Some hearing aid models do not support caller identification.

 Caller identification must be enabled by your hearing healthcare professional.

Using uStream with an audio device

Streaming music via Bluetooth

If you have an audio playback device equipped with Bluetooth functionality (e.g. MP3 player or computer), uStream may be used to stream audio wirelessly to your hearing aids.

Pairing your audio device

1. Enable Bluetooth on your device/phone.
2. Starting with the uStream in the off position, press and hold the (⊕) button and turn the uStream on. The blue audio indicator will flash.
3. Enable Bluetooth pairing mode on your device so that a new device can be added. See manufacturer instructions for details. If prompted, enter "0000" for pairing code.
4. After successful pairing, the blue audio indicator on the uStream will stop blinking.

Playing audio from your device

1. Begin playback on your device.
2. Adjust the volume level on the device to a medium to high level.
3. Make minor volume adjustments with the uStream volume buttons (+, -).

Pausing playback

If you want to temporarily stop streaming to hear a conversation, press the (A) button to pause the streaming. Press the (A) button again to resume streaming.

Stopping playback

To stop playback on your device and switch back to regular hearing programs, press and hold the access (A) button for two seconds.

 Note: Please refer to the device's user guide or dealer support if you cannot complete the pairing procedure with the generic steps described.


 Note: uStream streaming time is four hours with a fully charged battery.

 Note: During streaming, uStream cannot be used as a hearing aid remote control. Press and hold the access (A) button (two seconds) to stop the Bluetooth stream completely.

Using uStream with uTV 3

If you have a uTV 3, you can:

- Stream audio from your TV up to 30 m (90 ft.) to your hearing aids
- Recharge your uStream in the charging cradle

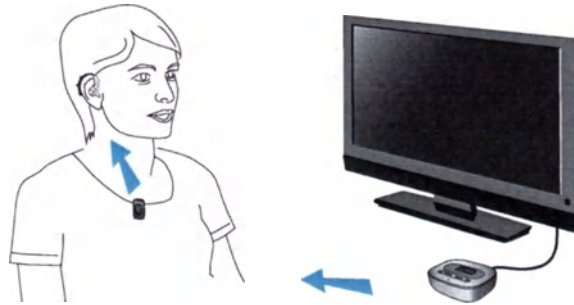Charging your uStream on the charging cradle

1. Ensure the uTV 3 has the correct charging cradle for the uStream – see uTV 3 manual for details.

2. Ensure the uTV 3 is turned on and place the uStream into the charging cradle.
3. Allow two hours for a full charge.

Pairing uTV 3 to uStream (not required for uTV 3 + uStream bundle)

1. Turn on the uTV 3 and press the pair button on the back of the uTV 3.
2. Starting with the uStream in the off position, press and hold the (+) button and turn the uStream on. The blue audio indicator will flash.
3. After successful pairing, the blue audio indicator on the uStream and the uTV 3 will stop blinking.

Playing back audio from your TV

1. Begin playback on your TV.
2. Adjust the volume level on the uTV 3 to a medium to high level.
3. Make minor volume adjustments with the uStream volume buttons (+/-).

Pausing playback

If you want to temporarily stop streaming to hear a conversation, press the (A) button to pause the streaming. Press the (A) button again to resume streaming.

Stopping playback


To stop playback on your device and switch back to regular hearing programs, press and hold the access (A) button for two seconds.

- ❗ Changing hearing aid programs is not possible while streaming.
- ❗ Please contact your hearing healthcare professional if you can not complete the pairing procedure with the steps described.
- ❗ To listen to stereo music from your device, it must support the A2DP/"Stereo Headset" profile.
- ❗ uStream is capable of streaming up to four hours with a fully charged battery.

Using uStream with uMic

If you have a uMic, you can:

- Stream audio from a person or device.

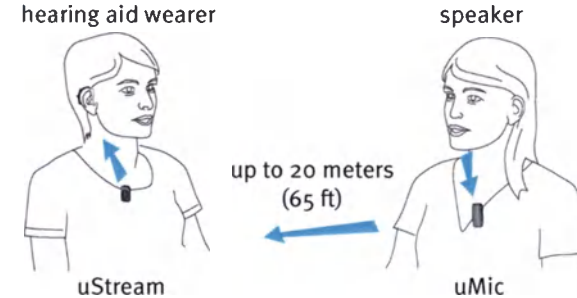Pairing uMic to uStream

1. Starting with the uMic in the off position, press and hold the (+) button and turn the uMic on. Continue to hold the (+) button until the blue audio indicator flashes.
2. Starting with the uStream in the off position, press and hold the (+) button and turn the uStream on. The blue audio indicator will flash.

3. After successful pairing, the blue audio indicator on the uStream and uMic will stop blinking.

Playing back audio from your uMic

1. Clip the uMic to a person, or position it near a sound source such as a speaker or television
2. Adjust the volume level on the uMic to a medium to high level
3. Make minor volume adjustments with the uStream volume buttons (+ -).

Pausing playback

If you want to temporarily stop streaming, press the (A) button to pause the streaming. Press the (A) button again to resume streaming.

Stopping playback

To stop playback on your device and switch back to regular hearing programs, press and hold the access (A) button for two seconds.

 Changing hearing aid programs is not possible while streaming.

 Please contact your hearing healthcare professional if you can not complete the pairing procedure with the steps described.

 uStream is capable of streaming up to four hours with a fully charged battery.

Managing multiple audio inputs to the uStream

When an audio source is connected via the uStream, your hearing aids will automatically switch to special streaming programs.

It is possible to have two Bluetooth audio sources connected.

Spoken status indicators and the indicator light inform you about the connected inputs.

- 📖 A phone call always has priority. You will hear the ringtone even when listening to other audio sources. If you reject or terminate the call, audio streaming will resume.
- 📖 It is not possible to change hearing aid programs using the access (A) button as long as any audio source is connected to the uStream.
- 📖 If you want to listen to music from your connected mobile phone, you must start the music via your phone initially; it cannot be started from the uStream.


Selecting from multiple inputs using the uStream

With only one audio source (e.g., from a Bluetooth transmitter), each press of the access (A) button will pause or resume the audio source. You can completely stop the music by pressing and holding the access (A) button for two seconds. It must be re-started using the music device.

Multiple sound sources can be accessed sequentially by pressing the access (A) button. After each press, the uStream selects the next connected source in a fixed sequence.

If a source is not connected, the uStream will skip over it and proceed to the next source.

Every connected audio source has a fixed place in the sequence, as shown in the diagram.

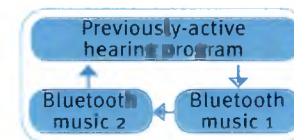Deleting the pairing of Bluetooth devices

If you encounter issues with one of the paired devices, or if a device stops connecting to your uStream, you may consider deleting the pairing of devices to the uStream.

1. Turn on your uStream while pressing the volume up (+) button.
2. Release the volume up (+) button only after ten seconds when the rapidly blue blinking indicator light turns green.
3. Also delete the uStream pairing entries on your mobile phone before any re-pairing.
4. Repeat the pairing procedure.

 **Deleting the pairing of Bluetooth devices** removes the ability to connect to these devices. You will need to perform the pairing process again if you wish to reconnect.

Troubleshooting guide

Cause	Possible remedy
During Bluetooth pairing my other device is prompting for a code, what do I have to enter?	
The pairing code is required for the initial coupling of uStream with your Bluetooth device (pairing)	The uStream pairing code is "0000" (four zeros)
I hear intermittency during audio streaming via uStream	
The Bluetooth device is too far away, or worn in a back pocket (body shielding)	Do not to exceed the maximum distance 5-10 m (15-30 ft) and keep Bluetooth devices in front of you
uStream searches for known devices for two min after switching on	After two minutes, uStream will stop searching and the intermittency will stop
A phone connected to other devices, such as a PC, may not be able to continuously stream music	When using your Bluetooth phone as music player, disconnect it from other devices, such as PCs

Cause	Possible remedy
The volume of the music or phone is too loud or too quiet	
Different audio sources have different outputs	To change the volume of the phone or any other sound source to uStream, first use the volume button of that external device, then make minor adjustments with the uStream
I have difficulties understanding the phone conversation in a noisy environment	
Your hearing aid microphones are set too high during the call	Reduce the environmental noise by pressing the volume down button for two seconds and then use the phone volume button to increase the phone volume
The caller can hear me but I cannot hear the caller	
The call was transferred to your phone	Make sure the call is transferred back to uStream via the phone menu or long press of the access (A) button
Some phones, when accepting the call through the phone button, don't use uStream	Always accept the call by pressing the access (A) button

Cause	Possible remedy
uStream is no longer recognized by the phone or other Bluetooth device to which it was previously paired	
Two minutes after no device was found, uStream stops seeking devices to save power	Switch uStream off, wait until the red indicator switched off, then turn device on again
The devices are too far away from each other	Bring devices within one meter of each other and try again
uStream can be paired to a maximum of six different Bluetooth devices. If the internal memory is full, new pairings may overwrite previous pairings	Repeat the pairing between uStream and the Bluetooth device which was over-written
Pairing was deleted	Repeat the pairing between uStream and the Bluetooth device which was deleted
I am not streaming but the hearing aids keeps switching between streaming and the normal hearing aid program	
Your phone sends sounds to uStream due to incoming SMS or mails	Disable the system sounds for reminders, alarms, SMS, etc. via your phone menu
Phone key tones are switched on	Switch off phone key tones via your phone menu

Cause	Possible remedy
I hear my phone's key-press tones in my hearing aids	
The key-press tones may be enabled on your phone	Deactivate key-press and confirmation tones on your phone. Please refer to the phone's user guide
Every time I switch on my phone, the music player starts	
For some phones this is normal behavior and not initiated by uStream	Some phones can be configured not to start the music player automatically. Press and hold the access (A) button to stop streaming. Or, stop the music player via phone menu after the phone is connected
When I get into my car, the uStream connection to my phone shows unexpected behavior	
The car Bluetooth system connects to your phone	Consider disconnecting the car system from your phone
I accepted the incoming phone call but cannot hear the call in my hearing aids	
Some phones don't use uStream when accepting the call through the phone button	Always accept the call by pressing the uStream access (A) button

Cause	Possible remedy
The phone is ringing but the indicator light is not solid blue and there is no ringing signal through the hearing aids	
Bluetooth is disabled in the phone	Ensure Bluetooth is enabled according to your phone's user guide
uStream and your phone are not connected to each other	1. Disable the Bluetooth functionality in your phone and turn it on again 2. Restart uStream by switching it off and then on. The uStream Bluetooth functionality will be automatically reactivated 3. Activate uStream by selecting it in the phone's "Active device" menu 4. Reduce the distance between uStream and the phone
uStream is not paired to the phone	Follow the Bluetooth pairing process in this guide
The distance is too great between uStream and the phone	Wear uStream clipped to your clothing. Keep distance to phone within <10m (<30ft)

Cause	Possible remedy
The phone is ringing, the indicator light is solid blue but there is no ringing signal through the hearing aids	
uStream is out of range of the hearing aids	Wear uStream correctly clipped to your clothing with maximum 40 cm (16") distance to the hearing aids
Phone is set to vibrate	Enable ring tone on phone
If you are already on a phone call, then a second call cannot be recognized	Hang up the first phone call, then accept the second call
My phone rings but the name of the caller is not being spoken by uStream	
The caller identification feature is not available or not enabled on your phone	Please refer to the phone's user guide to check the availability of this feature
Two phones are connected to your uStream	If two phones are connected, caller identification is only available on the phone that was last paired to uStream
Caller identification may be switched off	Please ask your hearing healthcare professional to enable caller identification on your uStream

Cause	Possible remedy
During a phone call the other party has difficulty hearing me	
uStream microphone opening may be covered	Make sure the microphone openings are not covered by part of your body, clothing or dirt and debris
uStream may not be worn correctly	Make sure your uStream is clipped securely to your clothing and within a 40 cm (16") distance to your hearing aids
uStream may rub against clothing	Reduce movement while you speak
The surrounding noise may be too loud	Bring uStream closer to your mouth. Suggest to your callers that they increase their phone volume
uStream is unresponsive or there is no indicator light when switching on	
uStream may be completely discharged	Charge your uStream for at least one hour
It may indicate a software problem	Switch uStream off, wait 30 seconds, then switch on again

Cause	Possible remedy
The operating time of your uStream is reduced substantially	
The battery typically needs several charging cycles to achieve full performance	Charge for at least three hours the first time you charge uStream. Let a new uStream discharge completely three times
I have interruptions while my audio streaming is active	
uStream is too far away from the hearing aids while streaming	Wear uStream clipped to your clothing with a maximum 40 cm (16") distance to the hearing aids
Bluetooth signal is out of range	Keep the Bluetooth device within the operating range
I cannot hear the demo sound	
Hearing aids are not worn correctly	Make sure your hearing aids are inserted correctly
uStream is not worn correctly or switched off	Wear uStream clipped to your clothing with a maximum 40 cm (16") distance to the hearing aids. Make sure uStream is charged and switched on
For any problems not listed in this user guide, please contact your hearing healthcare professional.	

Warnings

-  Please read the information on the following pages before using your Unitron uStream.

Hazard warnings

-  Keep this device out of reach of pets and children under 3 years.
-  Only use hearing aids that have been specially programmed for you by your hearing healthcare professional.
-  Changes or modifications to any of the devices that were not explicitly approved by Unitron are not permitted.
-  For persons with active implantable medical devices (i.e. pacemakers, defibrillators, etc.) caution should be exercised when using this device. Ask your hearing healthcare professional for details.
-  Opening the uStream might damage it. If problems occur which cannot be resolved by following the remedy guidelines in

the troubleshooting section of this user guide, consult your hearing healthcare professional.

- ⚠ This product has an embedded, non-replaceable battery. Do not attempt to open the product or remove the battery as this may cause injury and damage the product.
- ⚠ Dispose of electrical components in accordance with your local regulations.
- ⚠ Do not use the device in explosive areas (mines or industrial areas with danger of explosions, oxygen rich environments or areas where flammable anesthetics are handled) or where electronic equipment is prohibited.
- ⚠ External devices may only be connected if they have been tested in accordance with corresponding IECXXXXX standards. Only use accessories approved by Unitron.
- ⚠ Using your uStream cables in any way contradictory to their intended purpose (e.g., wearing the USB cable around the neck) can cause injury.

48

- ⚠ When operating machinery, ensure that no parts get caught in the machine.
- ⚠ The mini-USB port is to be used for the described purpose only.
- ⚠ Caution: electric shock. Do not insert plug alone into electrical outlets.

- ⚠ For safety reasons, only use chargers supplied by Unitron or stabilized chargers with a rating of 5VDC, min. 500 mA.

Information on product safety

- ⚠ Check your hearing aid compatibility with your hearing healthcare professional.
- ⚠ Protect the uStream from excessive moisture (bathing, swimming), heat (radiator, car dashboard), and direct skin contact when sweating (workout, fitness, sports).

49

- ⚠ Do not drop the device. Dropping onto a hard surface can damage your device.
- ⚠ Special medical or dental examination including radiation described below, may adversely affect the correct functioning of your remote. Remove and keep it outside the examination room/area before undergoing a medical or dental examination with X-ray (also CT scan) or any medical examinations with MRI/NMRI scans, generating magnetic fields.
- ⚠ Protect the device, charging slot and power supply from dirt and debris.
- ⚠ Clean the uStream using a damp cloth. Never use household cleaning products (washing powder, soap, etc.) or alcohol to clean the uStream.
- ⚠ Never use a microwave or other heating devices to dry the uStream.

- ⚠ The digitally-coded, inductive transmission technology used in the uStream is highly reliable and experiences virtually no interference from other devices. It should be noted, however, that when operating the hearing system near computer equipment, larger electronic installation or other strong electromagnetic fields, it may be necessary to be at least 60 cm (24") away from the interfering device to ensure proper operation.
- ⚠ Keep the device at least 10 cm away from any kind of magnets.
- ⚠ For safety reasons, recharge the uStream only with chargers supplied by Unitron or with USB-certified chargers >500 mA.
- ⚠ Do not connect an USB or audio cable exceeding 3 meters (9ft) in length to the uStream.
- ⚠ When the uStream is not in use, turn it off and store it safely.

Other important information

- ⚠ High-powered electronic equipment, larger electronic installations and metallic structures may impair and significantly reduce the operating range.
- ⚠ If the hearing aids do not respond to the uStream because of an unusual field disturbance, move away from the disturbing field.

Information and explanation of symbols

With the CE symbol, Unitron confirms that this Unitron product – including accessories – meets the requirements of the Medical Devices Directive 93/42/EEC as well as the Radio Equipment Directive 2014/53/EU on radio and telecommunications transmitters.

The numbers after the CE symbol correspond to the code of certified institutions that were consulted under the above-mentioned directives.

This symbol indicates that it is important for the user to read and take into account the relevant information in this user guide.

This symbol indicates that it is important for the user to pay attention to the relevant warning notices in this user guide.

Important information for handling and product safety.

Australia and New Zealand EMC and Radiocommunications compliance label.

Copyright symbol

This symbol shall be accompanied by the name and the address of the authorised representative in the European Community.

This symbol shall be accompanied by the name and the address of the manufacturer (who are placing this device on the market).

Operating conditions:

This device is designed such that it functions without problems or restrictions if used as intended, unless otherwise noted in these user guides.

0° to 45° Celsius and relative humidity of <95% (non condensing).

Transport conditions:

Temperature: -20° to +60° Celsius (-4° to +140° Fahrenheit).

Humidity: Up to 90% (non condensing).

The symbol with the crossed-out garbage bin is to make you aware that this device may not be thrown away as normal household waste. Please dispose of old or unused devices, at waste disposal sites intended for electronic waste, or give your device to your hearing healthcare professional for disposal. Proper disposal protects the environment and health.

The Bluetooth® word mark and logos are registered trademarks owned by Bluetooth SIG, Inc. and any use of such marks by Unitron is under license. Other trademarks and trade names are those of their respective owners.

Google Play is a trademark of Google Inc.

Apple and the Apple logo are trademarks of Apple Inc., registered in the U.S. and other countries. App Store is a service mark of Apple Inc.

Compliance information

Declaration of Conformity

Hereby, Unitron Hearing declares that this product meets the requirements of the Medical Device Directive 93/42/EEC and complies with Radio Equipment Directive 2014/53/EU. The full text of the Declarations of Conformity can be obtained from the manufacturer.

Notice 1:

This device complies with Part 15 of the FCC Rules and with RSS-210 of Industry Canada. Operation is subject to the following two conditions:

- 1) this device may not cause harmful interference, and
- 2) this device must accept any interference received, including interference that may cause undesired operation.

Notice 2:

Changes or modifications made to this device not expressly approved by Unitron may void the FCC authorization to operate this device.

Notice 3:

This device has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules and ICES-003 of Industry Canada. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This device generates, uses and can

radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this device does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

Reorient or relocate the receiving antenna.

Increase the separation between the device and receiver.

Connect the device into an outlet on a circuit different from that to which the receiver is connected.

Consult the dealer or an experienced radio/TV technician for help.

Australia: Supplier Code Number N15398

New Zealand: Supplier Code Number Z1285

Patient feedback

Record your specific needs or concerns and bring to your first office visit after getting your hearing aids.

This will help your hearing healthcare professional to address your needs.

☐

☐

☐

☐

☐

Additional notes

Unitron Hearing, a division of
National Hearing Services Inc.
20 Beasley Drive, P.O. Box 9017,
Kitchener, ON N2G 4X1 Canada

Sonova AG
Laubisrütistrasse 28
CH-8712 Stäfa, Switzerland

For a listing of Unitron group
companies, please visit
www.unitron.com

Appendix 5 – 029-5993-02c_uStream_Gd_ENG

Distributor

Unitron, Inc. 029-5993-02 © 2011-2017 Unitron. All rights reserved.

unitron.com

Appendix 6 – english-completeguide-icubeii

Amir Panah

Complete guide to iCube II

Wireless programming

unitron. Hearing matters

Appendix 6 – english-completeguide-icubeii






Table of contents

Quick reference page.....	3
Indicator lights	4
Overview	5
First time setup.....	6
Connecting iCube II to your computer.....	7
Battery status.....	8
Operating time.....	8
Using iCube II with a patient.....	8
Using iCube II without a patient.....	9
Important information.....	10
Troubleshooting guide.....	12
Declaration of conformity.....	13
Notices	13

Quick reference page

The following page can be removed from this booklet for use as a quick reference guide.

Indicator lights

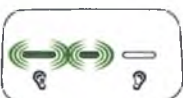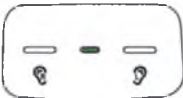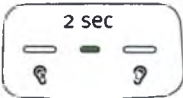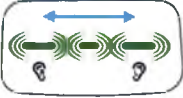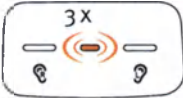Middle indicator remains solid orange		Charging	Right and middle indicator rapidly blinking for 2 seconds		Programming right hearing instrument
Middle indicator remains solid green		Charging complete	Left and middle indicator rapidly blinking for 2 seconds		Programming left hearing instrument
Middle indicator is solid green for two seconds		iCube II starting up	All indicators rapidly blinking back and forth		Pairing mode enabled
Middle indicator slowly blinking three times		iCube II switching off	All three indicators rapidly blinking		Updating firmware on iCube II
Middle indicator slowly blinking		On, no active connection to fitting PC	Middle indicator remains solid red		Malfunction, please refer to troubleshooting section
Middle indicator regularly blinking		Active connection to fitting PC			

Overview

iCube II enables the wireless transfer of fitting data from your computer directly to your patient's hearing instruments.

Your iCube II is compatible with Era™ and Unitron North platform hearing instruments.


For further information, please visit the Unitron website: www.unitron.com

Use the diagram to identify the components of the iCube II programmer.

iCube II accessories

- iCube II desktop stand
- USB hub

Contents may vary depending on country.

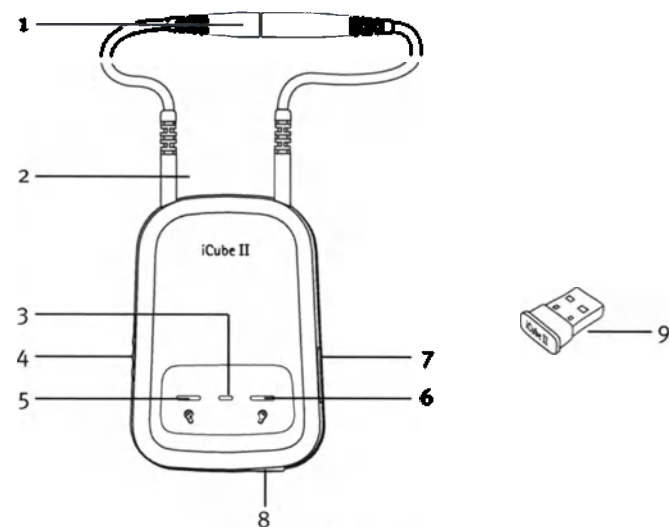Legend

- 1 Rear neckloop fastener
- 2 Neckloop plugs
- 3 Middle indicator (battery state and PC data transfer)
- 4 On/Off button
- 5 Right indicator light
- 6 Left indicator light
- 7 Pairing button
- 8 Mini USB charging input
- 9 Wireless iCube II USB adapter

First time setup

Charging the iCube II prior to setup


Connect the iCube II charger by plugging the power supply into a power outlet and the mini USB plug into the iCube II as shown. Alternatively you can charge your iCube II with the USB cable using the USB port on your computer system.

If you use the iCube II desktop stand, install it as shown.

Your iCube II battery can be charged hundreds of times. If the operating time of your iCube II reduces substantially, please contact your local Unitron representative.

Prior to the first use, ensure that the iCube II battery is fully charged.

Fully charged, your iCube II will operate for up to 12 hours of continuous use.

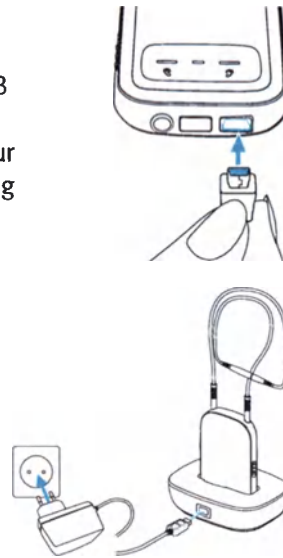 **Note:** When charging your iCube II for the first time, leave it to charge for at least three hours (even if the battery indicator starts flashing earlier).


 **Note:** Charging a depleted battery will take approximately two hours. It is safe to leave iCube II connected for longer periods, even if not in use.

 **Note:** A new battery's full performance is achieved only after three complete charge and discharge cycles.

Turning iCube II on/off

To turn on your iCube II, slide the switch to the on or off position, as shown. Once your iCube II is fully operational, the middle indicator will blink slowly.

Switch iCube II off when not in use and while charging. The middle indicator light will stop blinking.

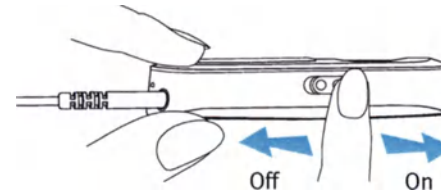Connecting iCube II to your computer

Your iCube II must be connected to your computer before it can be used for fitting.

There are three methods that can be used to connect the iCube II to your computer.

A. Using a wired USB connection

Using a USB cable, connect the iCube II to a spare USB slot on your computer. This is recommended when performing technology or firmware upgrades.

B. Using the USB adapter

Plug the pre-paired USB Bluetooth® adapter into a spare USB port on your computer.

Put the iCube II into pairing mode by pressing the  button on the side until the three lights flash.

C. Using the Bluetooth connection in your computer

To use iCube II wirelessly via Bluetooth, make sure that you have the Bluetooth function enabled on your computer. For specific instructions, please refer to the

documentation of the computer and operating system that you are using.

If you are using NoahLink and you have already installed the corresponding external Bluetooth adaptor, you can use that same adaptor for iCube II.

Refer to the “quick start guide” document available with the adaptor’s software CD and follow the instructions.

 **Note:** iCube II works within a Bluetooth operating range of up to 10 meters (33 ft.) of the fitting computer.

Put the iCube II into pairing mode by pressing the pairing button  on the side until the three LED lights flash.

The final step is to enable your iCube II in Unitron TrueFit™ fitting software. This is done by enabling iCube II support in [Options > Preferences > Fitting Device > iCube](#) and selecting a connection type.

Detailed iCube II instructions for setting up Unitron TrueFit with iCube II can also be found in the guides section of Unitron TrueFit help.

Removing a Bluetooth pairing

If you wish to delete the Bluetooth pairing from your iCube II, press and hold the pairing button for about ten seconds until the middle indicator is lit solid blue.

 **Note:** When connecting iCube II for the first time via USB, your computer will start searching for the associated driver software. Depending on the operating system, you may be presented with an installation wizard: Found New Hardware Wizard. If so, let your computer automatically install the necessary driver software. The driver software has already been installed during the fitting software installation and should be identified by the wizard. It is not necessary to insert the installation CD.

Battery status

While turned on, your iCube II indicates the battery status on the middle indicator as follows:

- Solid green – Battery level is high
- Solid orange – Battery level is low

For safety reasons, recharge the battery only with chargers supplied by Unitron.

Normally, charging an empty battery takes about two hours. It is safe to leave the charger connected overnight – the battery cannot be overcharged.

 **Note:** A new battery's full performance is achieved only after three complete charge and discharge cycles.

 **Note:** The iCube II battery can be charged several hundred times. If the operating time of your iCube II reduces substantially, please contact your local Unitron representative.

Operating time

Fully charged, your iCube II will operate for up to 12 hours of continuous use.

Using iCube II with a patient

Batteries need to be inserted into the hearing instruments and the hearing aids need to be turned on prior to detection.

Before using iCube II to fit a patient with hearing instruments, make sure the battery is sufficiently charged. Unplug the charging cable from the iCube II and turn it on.


Make sure that iCube II is selected and shown as the fitting device in Unitron True:fit.

Wearing iCube II

Make sure both ends of the neckloop are firmly plugged into the iCube II.


Open the neckloop by unplugging the connector in the middle of the neckloop and place iCube II around your patient's neck.

Close the neckloop to allow the built-in antenna to operate correctly.

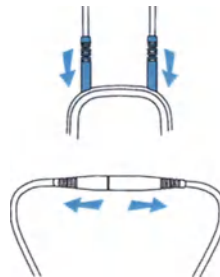Using iCube II without a patient

When programming hearing instruments without having your patients wear them, make sure they are turned on, and are in range of the iCube II. Please note that they are best positioned inside the loop.

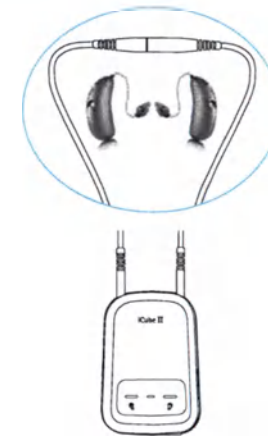Important information

Maintenance and care

- ⚠ Clean iCube II using a damp cloth. Never use household cleaning products or alcohol to clean the device.
- ⚠ Protect iCube II from excessive moisture, heat, excessive shock and vibration.
- ⚠ When iCube II is not in use, turn it off and recharge it so it will be sufficiently charged and ready for the next fitting session.
- ⚠ Do not drop the device. Dropping it onto a hard surface can damage the device.
- ⚠ Do not disconnect the neckloop while the device is transmitting signals to the hearing aids.

Safety notice and other important information

- ⚠ Keep this device out of reach of children and pets.
- ⚠ Persons using medical devices such as pacemakers are not permitted to use this device.
- ⚠ X-ray radiation (e.g. CT scans, MRI scans) may adversely affect the correct functioning of this device. We recommend that you keep iCube II outside the X-ray radiation procedure rooms.
- ⚠ High-powered electronic equipment, larger

electronic installations and metallic structures may impair and significantly reduce the operating range.

- ⚠ External devices may only be connected if they have been tested in accordance with applicable IEC Standards.

Important information

- ⚠ The digitally-coded, inductive transmission technology used in this device is extremely reliable and experiences virtually no interference from other devices. It should be noted, however, that when operating the device near a computer terminal or other strong electromagnetic fields, it may be necessary to be at least 60 centimeters (24”) away to ensure proper operation.
- ⚠ If the hearing instruments do not respond to the device because of an unusual field disturbance, move away from the disturbing field.
- ⚠ Use only original Unitron accessories.
- ⚠ Opening iCube II might damage it. If problems occur which cannot be resolved by following the remedy guidelines in the troubleshooting section of this user guide, contact your local Unitron representative.
- ⚠ iCube II may only be repaired by an authorized service center.

Appendix 6 – english-completeguide-icubeii

- ⚠ Changes or modifications to the device that were not explicitly approved by Unitron are not permitted. Unitron is not responsible for situations resulting from a modified device.
- ⚠ If you use several iCube II devices in the same building, be aware that their inductive field can disturb one another if within a range of 2 meters (6 feet). This may even occur when used in different rooms.
- ⚠ Do not use iCube II while charging.
- ⚠ Caution: electric shock. Do not insert plug alone into electrical outlets.
- ⚠ This product has an embedded, non-replaceable battery. Do not attempt to open the product or remove the battery as this may cause injury and damage the device.
- ⚠ Do not use the device in explosive areas (mines or industrial areas with danger of explosions, oxygen rich environments or areas where flammable anesthetics are handled) or where electronic equipment is prohibited.
- ⚠ Dispose of electrical components in accordance with your local regulations.
- ⚠ Keep the device of at least 10 cm away from any kind of magnets.
- ⚠ Using the device, its connectors or cables in any way contradictory to their intended purpose as described in this user guide (e.g. wearing the USB cable around the neck) can cause injury.
- ⚠ Do not connect an USB cable exceeding 3 meters (9 feet) in length to the device.
- ⚠ When operating machinery, ensure that no parts get caught in the machine.
- ⚠ When the device is not in use, turn it off, keep the neckloop closed and store it safely.
- ⚠ The mini-USB port is to be used for the described purpose only.
- ⚠ High-powered electronic equipment, larger electronic installations and metallic structures may impair and significantly reduce the operating range.
- ⚠ For safety reasons, only use chargers supplied by Unitron or stabilized chargers with a rating of 5VDC, min. 500 mA.
- ⚠ Protect the device connectors, plugs, charging slot and power supply from dirt and debris. Never use a microwave or other heating devices to dry the device.

Appendix 6 – english-completeguide-icubeii

Troubleshooting guide

Cause	Possible remedy
iCube II is no longer recognized by the computer to which it was previously paired	
The iCube II USB adapter is disconnected from the fitting computer	Plug in the iCube II USB adapter
iCube II is out of range	Move iCube II closer to your computer
iCube II turns off	
The battery is empty	Charge iCube II
iCube II is unresponsive and the middle indicator is lit solid red	
This may indicate a malfunction	Switch iCube II off and wait for 20 seconds before switching on again. If iCube II is connected to the fitting PC via cable, disconnect it from the cable and perform reboot as described above

Cause	Possible remedy
A second or different iCube II cannot be connected to the fitting PC	
Only one iCube II can be paired with the fitting PC	If you use the USB adapter, press the pairing button on iCube II to connect it with your computer. If you do not use the USB adapter, make sure your Bluetooth device is switched on and in pairing mode to establish the connection

Declaration of conformity

Unitron Hearing declares under its sole responsibility that the iCube II is in conformity with the following standards or other standardizing documents:

Radio: EN 300-330
EN 300 328

EMC: EN 60601-1-2

Safety: EN/IEC 60601-1

Implementing the provisions of the MDD 93/42/EEC (Medical Device Directive), R&TTE Directive 1999/5/EC (Radio and Telecommunications Terminal Equipment), EMC Directive 2004/108/EC (Electromagnetic Compatibility) and LVD 2006/95/EC (Low Voltage Directive).

Notices

Notice 1

This instrument is certified under:

 FCC ID: KWC-ICube2
IC: 2262A-ICube2

This device complies with Part 15 of the FCC Rules and with RSS-210 of Industry Canada. Operation is subject to the following two conditions:

1. this device may not cause harmful interference, and
2. this device must accept any interference received, including interference that may cause undesired operation.

Changes or modifications made to this equipment not expressly approved by Unitron Hearing may void the FCC authorization to operate this equipment.

Notice 2

This Class B digital apparatus complies with Canadian ICES-003.

Appendix 6 – english-completeguide-icubeii

Notice 3

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules.

These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses, and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation.

If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna
- Increase the separation between the equipment and receiver
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected
- Consult the dealer or an experienced radio/TV technician for help.

The CE symbol is confirmation by Unitron that this Unitron product meets the requirements of directive 93/42/EEC on medical devices and the R&TTE Directive 1999/5/EC on Radio and Telecommunication Terminal Equipment.

The symbol with the crossed out disposal bin indicates that this product shall not be treated as household waste. Please hand over your old or unused product to the applicable collection point for the recycling of electrical and electronic equipment or bring your old product to your hearing care professional for appropriate disposal. By ensuring this product is disposed of correctly, you will help prevent potential negative effects on the environment and human health.

Operating conditions:

This product has been designed for trouble-free operation without restrictions when used as intended, unless otherwise stated in this user guide.

Transportation and storage conditions:

Temperature should not exceed limits of $-20^{\circ}/60^{\circ}$ Celsius at a relative air humidity of 65% for extended periods during transportation and storage. Air pressure between 500 and 1100 hPa is not detrimental to this device.

Appendix 6 – english-completeguide-icubeii

 Z099 Australia & New Zealand

 Brazil

 CMI ID: China

 Japan

 MIC South Korea

Appendix 6 – english-completeguide-icubeii

At Unitron, we care deeply about people with hearing loss.
We work closely with hearing healthcare professionals to provide
hearing solutions that improve lives in meaningful ways.
Because hearing matters.

unitron.com

Unitron Hearing
20 Beasley Drive, P.O. Box 9017,
Kitchener, ON N2G 4X1 Canada

Unitron Hearing GmbH
Max-Eyth-Straße 20, 70736
Fellbach-Oeffingen, Germany

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PROVINCE OF ONTARIO

**Regional Municipality
of Waterloo**

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
DO CERTIFY AND ATTEST that the paper-writing hereto annexed is a document purporting to be

a Unitron report written in Cyrillic script, approved and signed on May 16, 2018 by Amir Elahi-Panah,
Manager, Quality Assurance,

an act whereof being requested I have granted under my Notarial Form and Seal of Office to serve and
avail as occasion shall or may require.

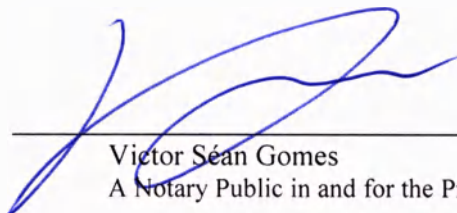IN TESTIMONY WHEREOF I have hereto subscribed my name and affixed my Notarial Seal of
Office at Kitchener, Ontario, in the Province of Ontario, this 23rd day of May, 2018.

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A/16-034 029-5905-26 ©2015 Unitron. All rights reserved.
A/16-034 029-5989-26 ©2015 Unitron. All rights reserved.
A/16-034 029-5990-26 ©2015 Unitron. All rights reserved.
A/16-034 029-5992-26 © 2015, Unitron. All rights reserved.
A/16-034 029-5993-26 ©2015 Unitron. All rights reserved.