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
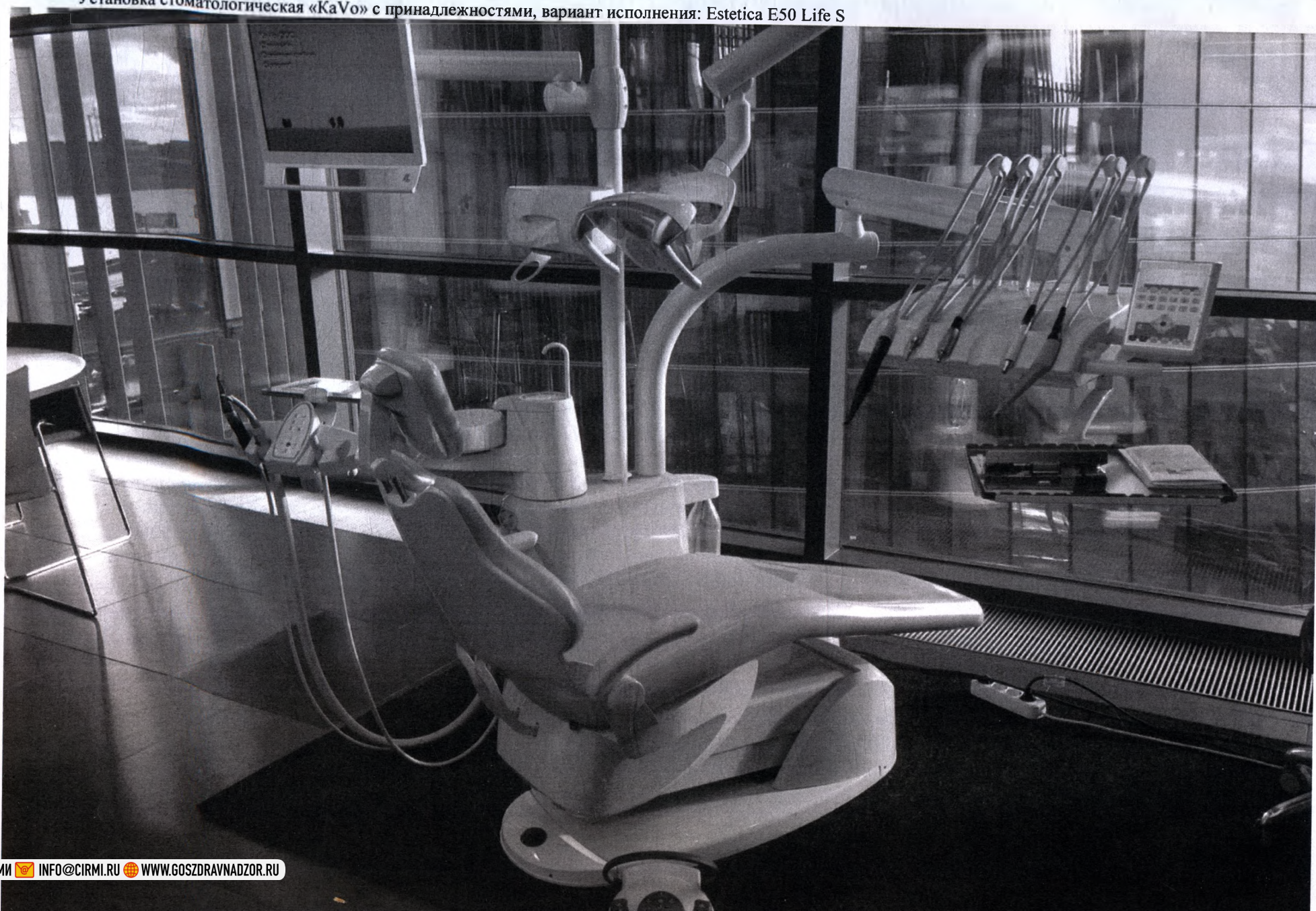
CONEXIO

## Installation Instructions/ Installationsanweisung


Insert the CONEXIO USB-Card into your PC to start the installation.

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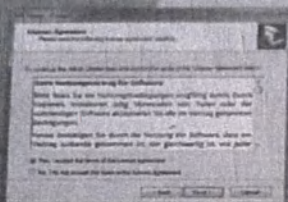2. Accept the license agreement to continue with the installation.

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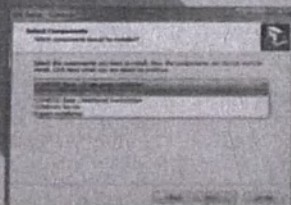3. Select the required installation type.

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
### Information


If a multi-site installation of CONEXIO (workplaces) is planned, firstly you have to install a work station with SQL-data base to connect each of the work stations with the common data base.  
Is it a single work station installation skip the step 4.

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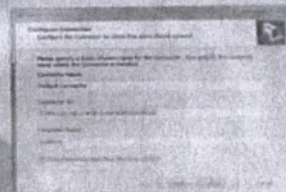
4. Is it an Auto-installation or manual installation? Enter the computer name of the SQL-Server to connect the Client with the Server Database.

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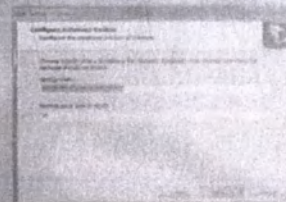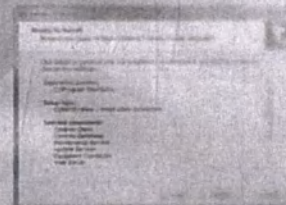
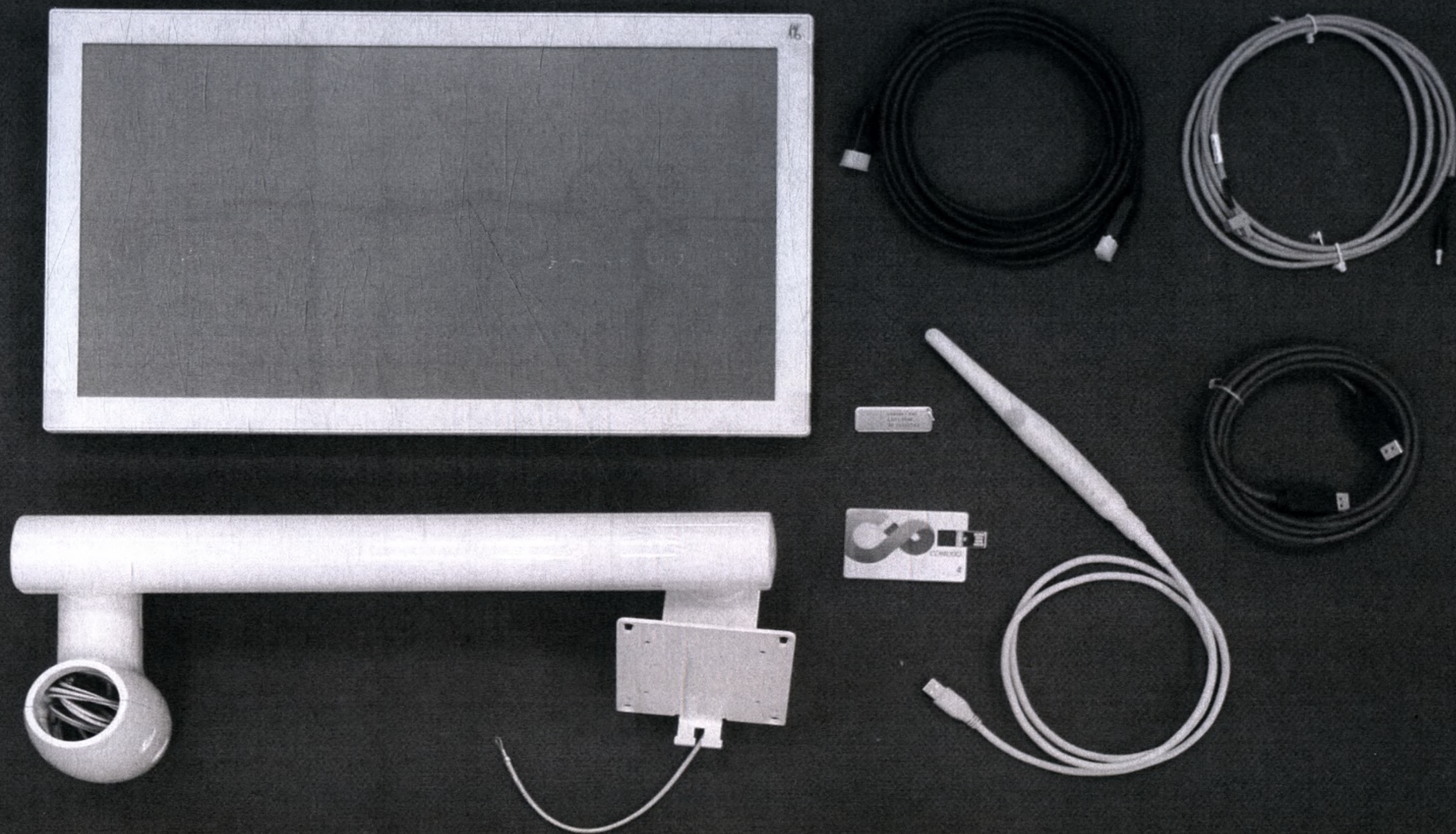
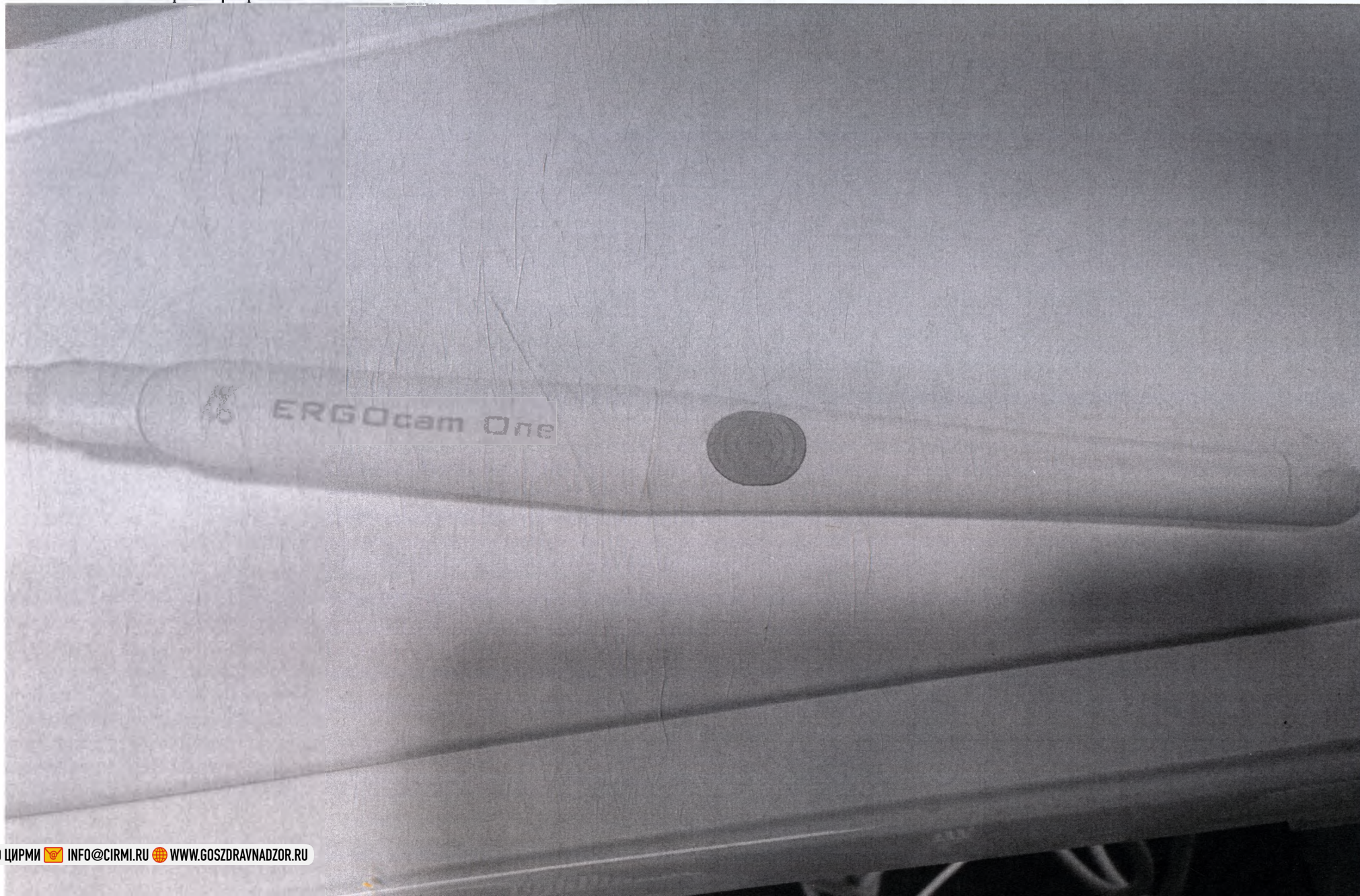
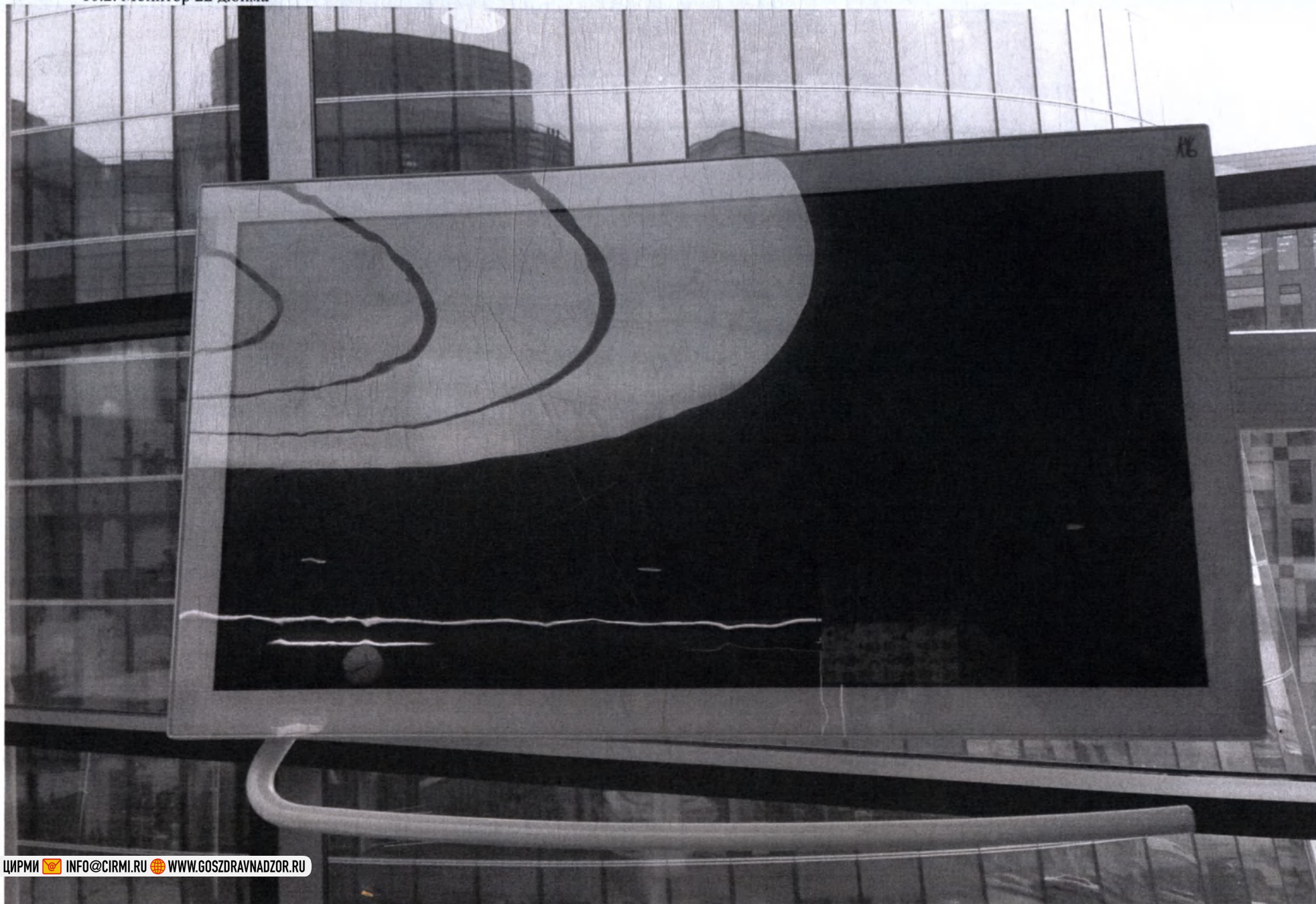
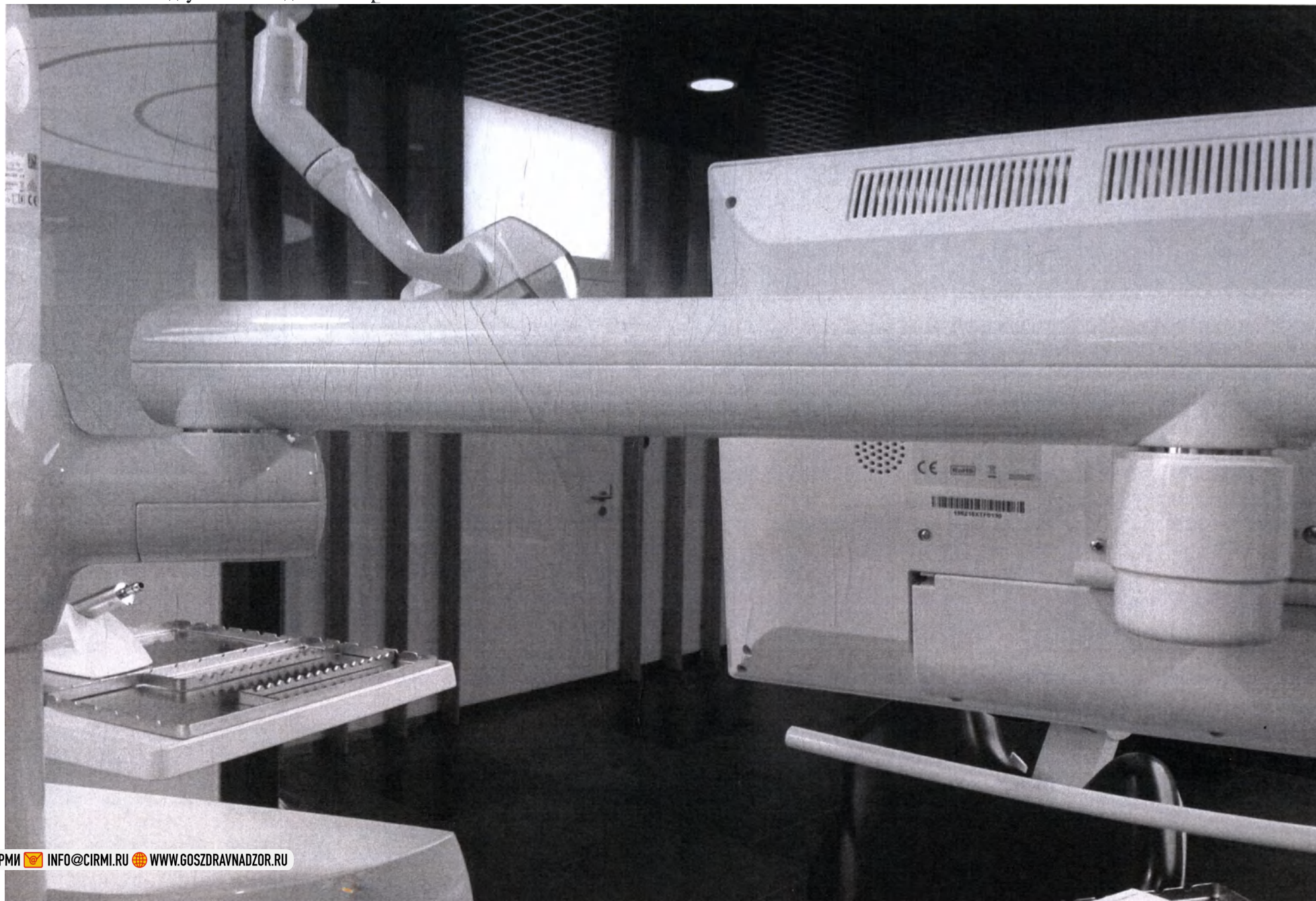
CONEXIO

## Installation Instructions/ Installationsanweisung


Insert the CONEXIO USB-Card into your PC, to start the installation.

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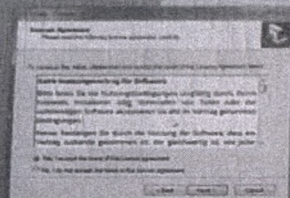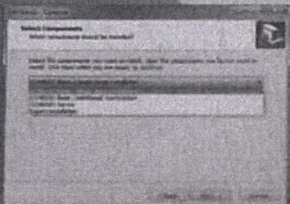
### Information


If a multi-user installation of DCHL XIS (and Manual work stations) is planned, firstly you have to install a server station with SQL-database to connect each of the work stations with the common data base. If it is a single work station installation skip the step 4.

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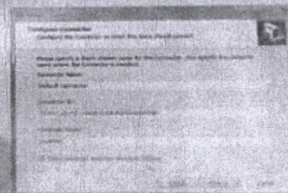
4. Is it an additional server station installation, when the connection to the computer Name of the SQL-Server is possible, the Client with the Server Database.

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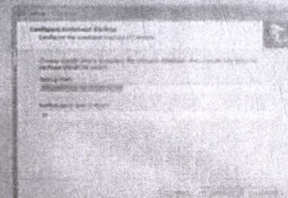5. Choose a required folder for the database location and select the backup folder.

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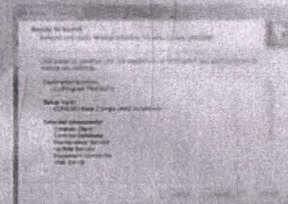
6. Choose a required folder for the backup location and select the backup folder.

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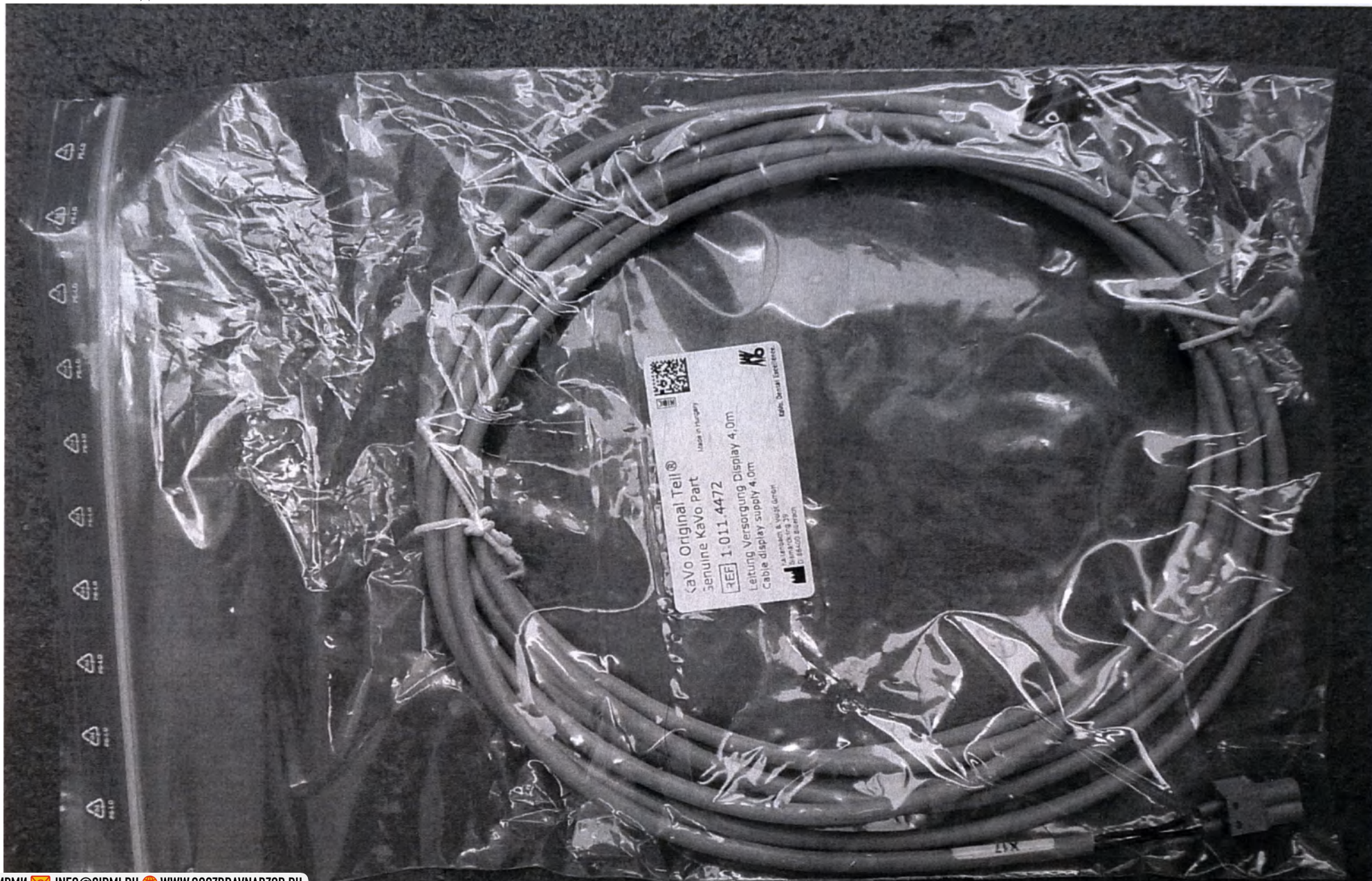
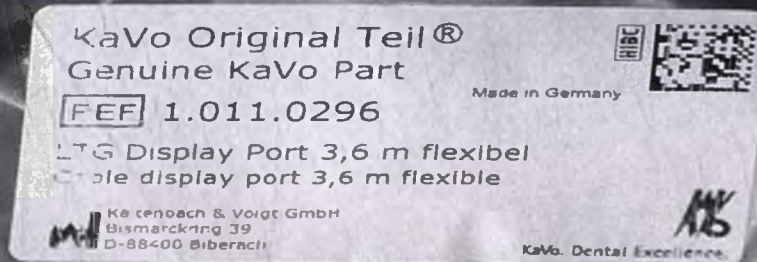
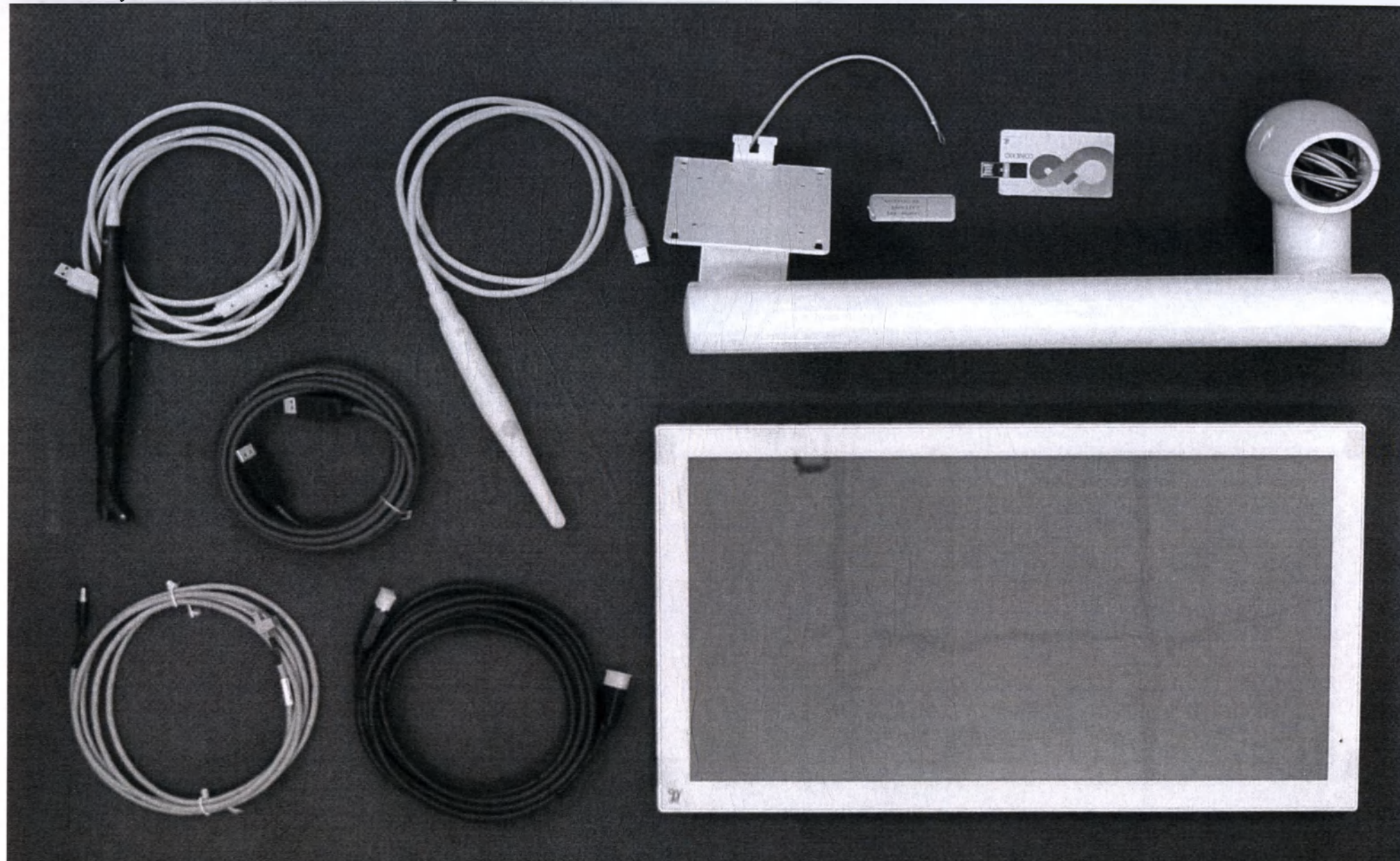
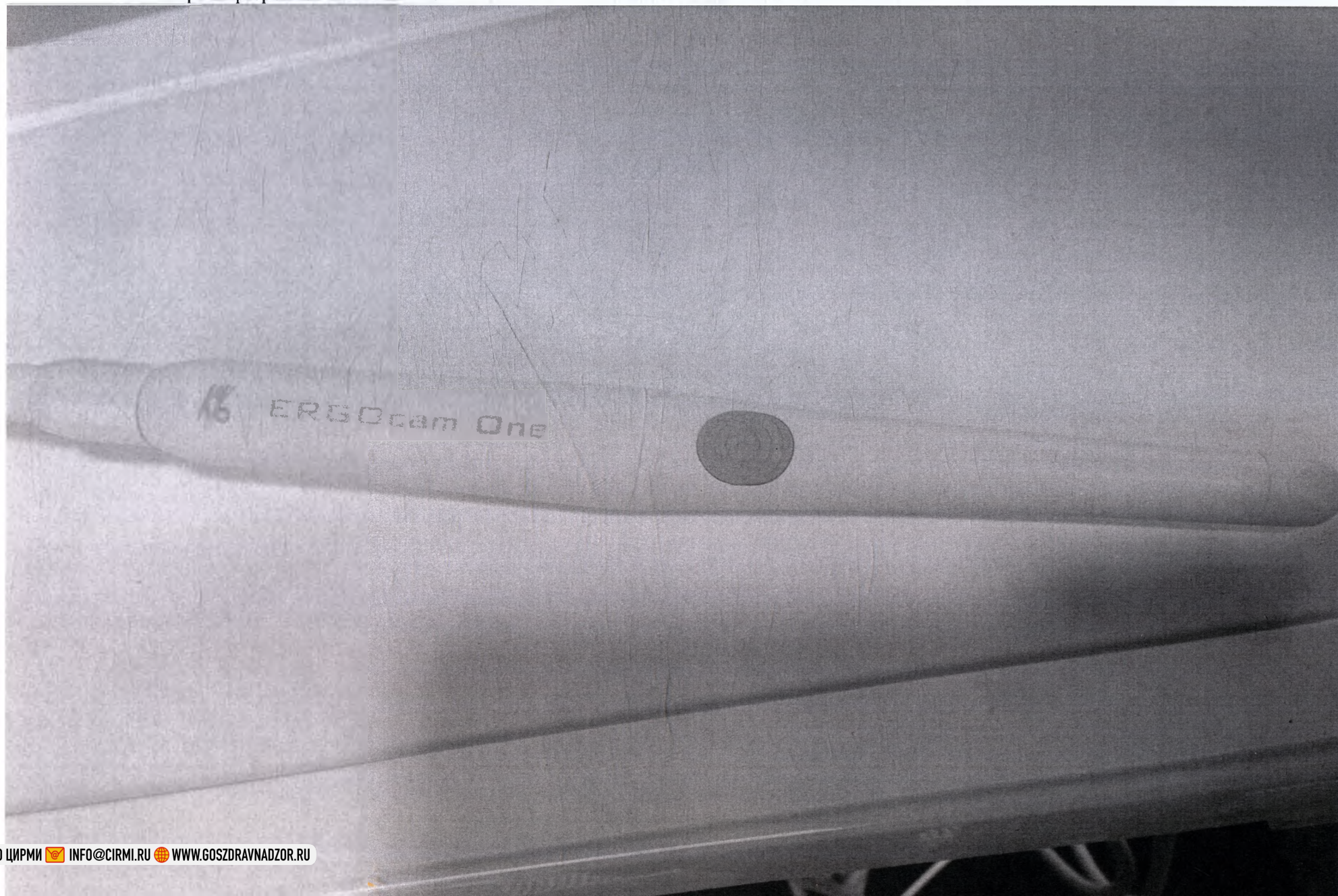
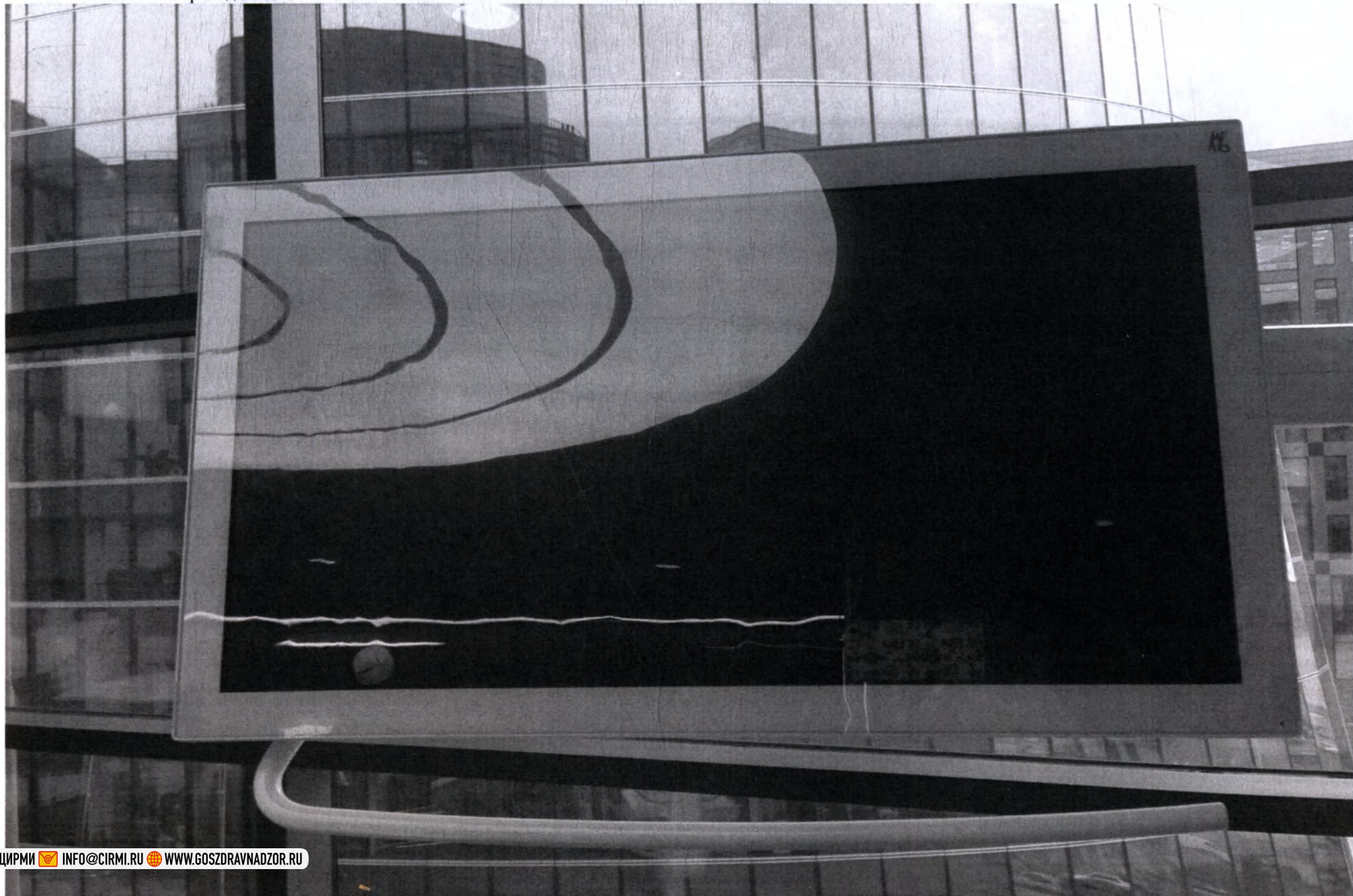
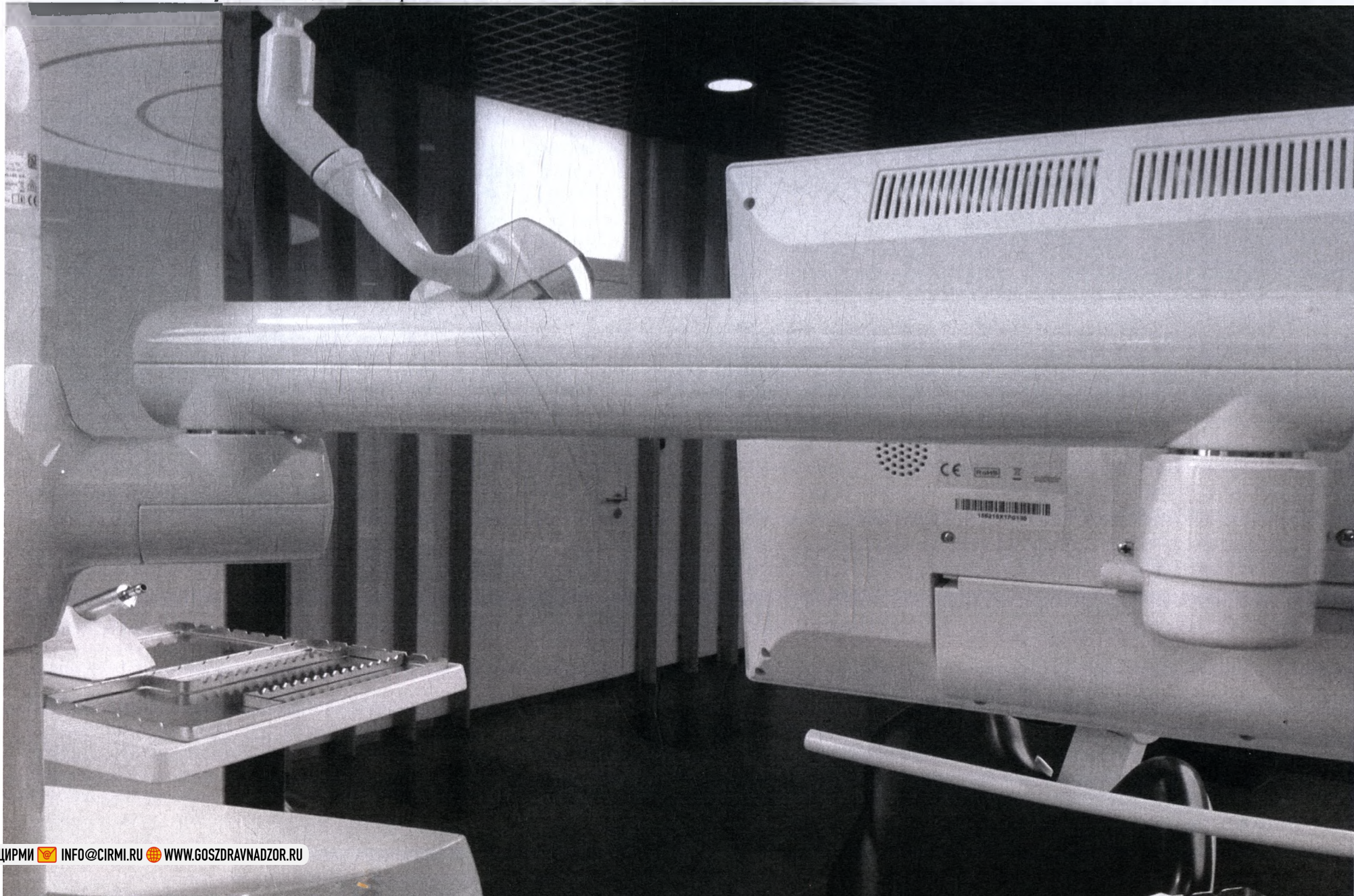
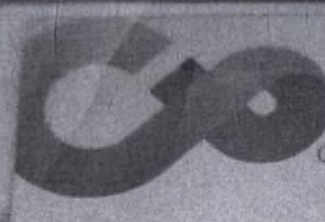
CONEXIO

## Installation Instructions/ Installationsanweisung


Connect the CONEXIO USB-Card  
to your PC to start the installation.

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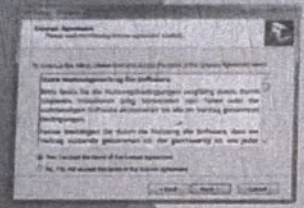2. Accept the license agreement to  
continue with the installation.

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
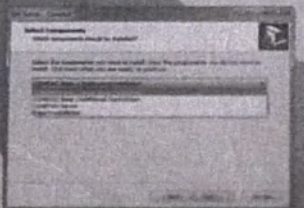
**Information:**  
If a multi user installation (Zusätzliche Arbeitsplätze) is planned, firstly you have to install a work station with  
SQL-data base to connect each of the work stations with the common data base.  
Is it a single work station installation skip the step 4.

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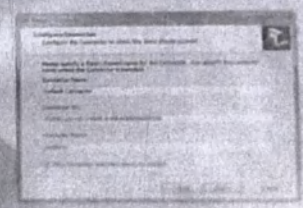
4. Is it an additional server installation, enter the computer name of the SQL-Server to connect the Client with the Server Database.

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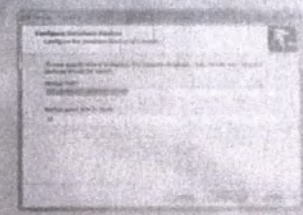5. Check a required folder for the  
database files, we select the  
default location.

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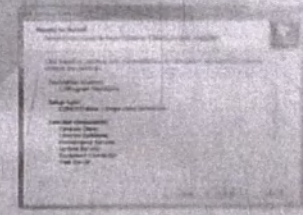6. Select a user who has the  
permissions to create the  
database in the server.

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