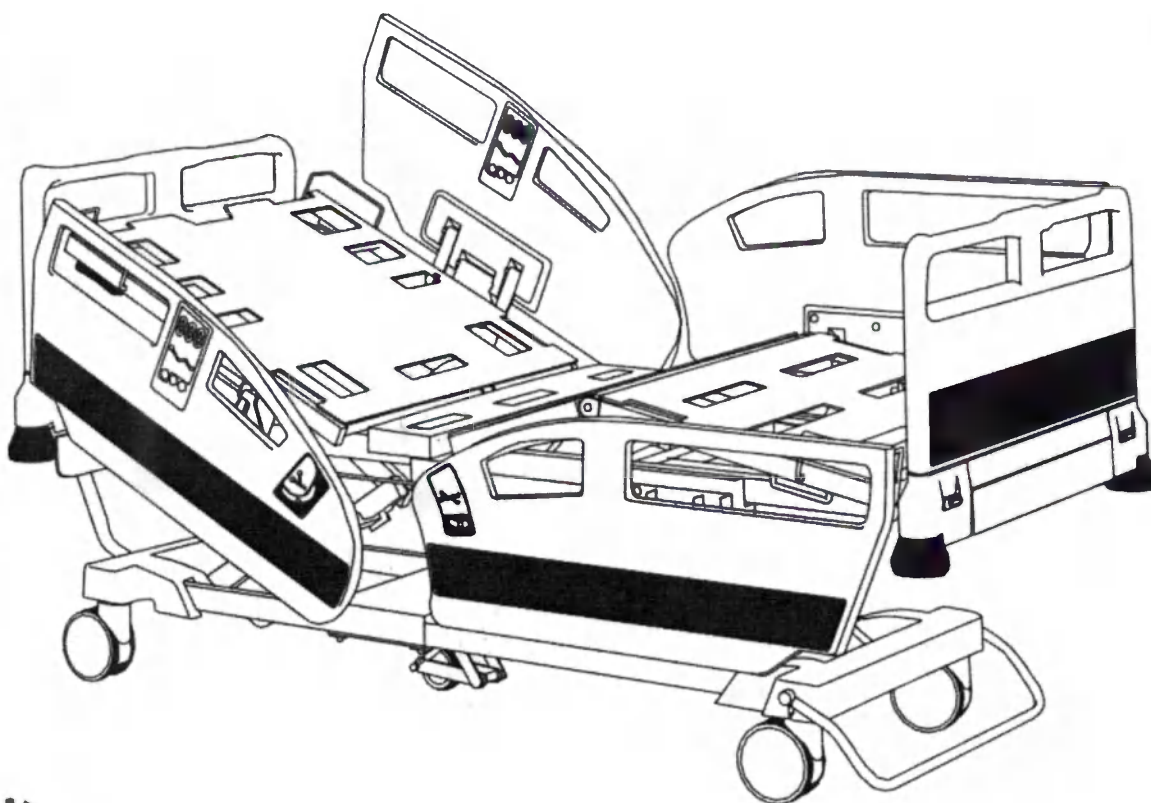
Medical functional electric 4-section NeWCare bed with accessories, versions NewCare Standard, NewCare S, NewCare XS, NewCare XSS

I

May2020| Cod XX8499900000522_ManualUso_NeWCare_EN_2020

User and Maintenance manual of NewCare beds.
Hospital bed

This manual is valid for the following bed models:

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www.pardo.es

Tel: 976 300 033 - Fax: 976 320 647

The manufacturer reserves the right to modify the content of this manual related on the product technical regulations.
For this reason, the manual content can have discrepancies with the real product

Index

SYMBOLS AND LABELS.....	5
Typographical conventions.....	6
Symbology.....	6
Organization.....	6
Bed labels.....	7
SECURITY ADVICES.....	8
General warnings.....	9
Electrical risk warnings.....	10
Entrapment warnings.....	10
REGULATIONS, RULES AND CODES.....	11
DESCRIPTION.....	13
Bed description.....	14
Bed components.....	16
Electrical components.....	16
Controls.....	16
Control box.....	16
Battery (optional).....	16
Linear electric actuators.....	17
Mechanical components.....	18
Side rails.....	18
Bed extension.....	18
Feet section position.....	18
Headboards and footboards.....	19
Wheels.....	19
Braking of bed.....	19
5 th wheel (optional):.....	19
RCP bilateral actioning system.....	20
Accessories support.....	20
Guides for X-ray cassette carrier.....	20
Protection discs.....	20
Accessories (optional):.....	20
IV-Pole.....	20
Trapeze bar.....	20
Balkan frame.....	20
Table under footboard.....	21
Multifunction table.....	21
Tank holder.....	21
Rx cassette tray.....	21
Light under bed.....	21
Brake alarm.....	21
List of components.....	22
BED OPERATION.....	24
Reception.....	25
Installation.....	25
Initial operation.....	25
Controls functioning.....	26
Kind of controls.....	26

Wired hand control:	26
Nursing controls:	27
Membrane controls:	28
Lifting pedal:	28
Operating	29
Basic control operating	29
Pre-programmed functions	30
Light indicators	31
Bed function locking	31
Other control functions	32
Bed components functioning	33
CPR bilateral actioning system	33
Rails	33
Full length rails	33
Tube side rails	34
Removable Headboard / Footboard with locking system	34
Placement	34
Removal	35
Multifunction pedal	36
Braking of bed	36
Moving of bed	36
Directional or 5 th wheel actioning	36
5 th wheel	36
Bed extension	37
Extension of bed	37
Removing from bed	37
Feet module positioning system	37
Extended table under footboard	38
Table extraction	38
Table folding	38
Cassette carrier	38
Bed base panels	39
Movement indicator	39
Backrest	39
Trendelenburg and reverse trendelenburg	39
CLEANING AND MAINTENANCE	40
Bed cleaning	41
Bed maintenance	42
Preventive maintenance guide	43
Rapid problems resolution	46
Product disposal	47
Warranty	47
TECHNICAL SPECIFICATIONS	48
Physical bed description	49
Rails physical description	50
Recommended mattress description	50
Use environment (according to direction EN 60601-2-52:2010)	50
Use environment physical description	50
Storage environment physical description	51
Electric system description	51
Electrical diagram	51
"EC" DECLARATION OF CONFORMITY	53

Anti-decubitus mattresses.

Technical characteristics

Variants.....

Density 100

Density 150

Density 250 fresh

Density 400

Density 450

Density 475 fresh

Density 600

Density 650 fresh

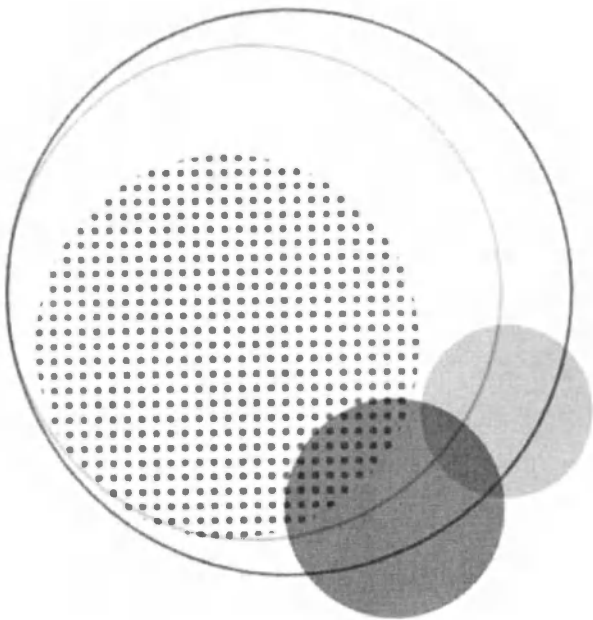
Density 900

Gigienic covers for mattresses.....

Cugar

Select

Matress clining

Symbols and labels

Symbols, labels and manual organization

Typographical conventions

This manual contains different typographical conventions and symbols designed to enhance readability and understanding of its content, for example:

- **Standard text:** used for the text throughout the manual.
- **Boldface text:** used to emphasize words or important phrases.

Symbology

- **Warning or general caution:**

It identifies situations or actions that may affect patient or user safety, as well as caution points out special procedures or precautions that the service personnel must follow to avoid equipment damage. Ignoring a warning may cause injuries to patient or user.

- **Electrical shock warning:**

Improper handling of some elements may cause electrocution by electrical exposure. Thus, it is recommended to be handled the marked areas by specialised personnel.

- **Entrapment and/or crushing warning:**

Do not hinder certain areas whilst moving as may cause entrapment or crushing.

- **Feet crushing warning:**

Do not hinder certain areas whilst moving to avoid crunching.

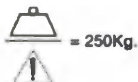
Organization

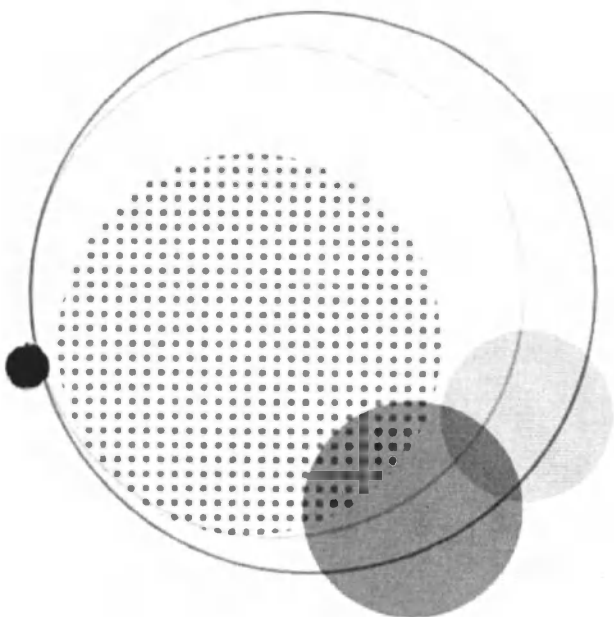
This manual has 3 chapters that includes from the most general issues to the most particular ones.

- **General description.**
Provides a general description of the product to know its main aspects.
- **Functioning and handling.**
Describes the use of different mechanical and electrical systems mounted in this product.
- **General procedures, cleaning and maintenance.**
Describes the cleaning tasks and preventive maintenance that have to be taken for a correct product running.

Bed labels

	Read the user manual
Input: (120V-220V) - 50/60Hz 700W	Motor Voltage
	Electric insulation Class II
	Protection against electric hazard accidents - Type B
	IPX6 Protection
	Suitable for indoor use only
	Safe working load (SWL)
	Maximum patient weight (MPW)
	Do not load over 140 kg at the beds ends whilst it is extended
	Serial number and CE mark label.
	Feet crushing warning
	Do not store any object under the bed
	Entrapment danger warning in these areas whilst lowering the bed
	Manual emergency CPR
	Recommended mattress dimensions
	Do not throw away the device along with the domestic waste. Recyclable product
	Do not seat on the articulated bed base panels whilst are raised, neither activate movements if there are people seated on the modules

 	Compliance to EU requirements of medical device
Symbols on the bed and package	
	Catalogue number
	Serial number
	Temperature range
Note : some symbols are not apply to all parts of device	
	Do not transport on a cart
	Keep away from moisture
	Industrias H. Pardo S.L. logo
	Fragile, handle with careБеречь от влаги
	Stacking is prohibited

Security advices

Security advices

General warnings

- This bed is not designed to be used with special anatomy users (amputees, etc) neither with children, because his anatomy or special circumstances can cause entrapment risks not contemplated in our risks tests. This is the reason why this product is not recommended for such cases.
- This manual describes all beds functions as well as its use, explaining all functions and components in detail. To obtain more accurate information relative to the bed handling refer to the manual specific sections of the manual.
- Refer to the manual before using the bed. Otherwise personal injuries or material damages may occur.
- The determination and problem solutions of the bed can be only done by qualified and authorised technical staff. The determination and problem solutions of the bed done by non authorised personnel can lead to personal injuries or material damages.
- Do not remove RFID labels from the bed. The reference that appears on it are essential to identified the bed and warranties the safety maintenance operations. Otherwise personal injuries or material damages may occur.
- Never raise the bed with other tool than the designed by the manufacturer. Not complying with this instruction may cause a rising mechanism dysfunction and body injuries.
- The motorised mechanisms of the machine can cause injuries. Do not use the bed controls until all the people are away from it. A bed articulation function can be interrupted releasing the control, pressing the button of the contrary movement or unplugging the supply connector from mains.
- The bed is designed for indoor use. Do not use or storage it on outdoors or humid environments.
- Do not use the equipment if the width, height and patient weight exceed the limits specified in this manual. Patient injuries may occur.
- Do not use in this bed mattress that cannot be adapted to the bed base articulations.
- Only authorised maintenance staff should make the preventive maintenance of the bed. A preventive maintenance mad by unauthorised personnel could produce personal injuries or material damages.
- The motors use cycles (intermittent operation) does not exceed 2 minutes / 20 minutes, otherwise equipment damages can be caused.
- Refer to the beds user manual before using the accessories to avoid any patient or user injuries risk.
- Only the following accessories specifically identified can be used in this bed. The use of any other accessories not included in the list will compromised the bed and patient safety.
- Do not hit the different movement controls. In the same way, never manipulate the inside of the manual or motor controls.
- Do not activate any movement before placing the mattress previously.
- Unplug the bed before any cleaning operation, bed moving or electric equipment repair is taking place. There is electrical shock risk.
- Never manipulate the manual or motor control inside.
- The electrical installation which the bed is connected has to have the regulatory protections, otherwise do not connect the device.
- In case any deformation is observed in any of the bed element, get in contact with the distributor.
- In the modules rising and lowering movements, it has to be observed that no materials or body parts are in-between the bed base and the bed supporting frames. Otherwise there are entrapment risks between the bed base and its supporting points whilst lowering.

- For the patient safety, the bed should be in its lowest position when the patient is not being attended, to avoid the risk of falling damages.

Electrical risk warnings

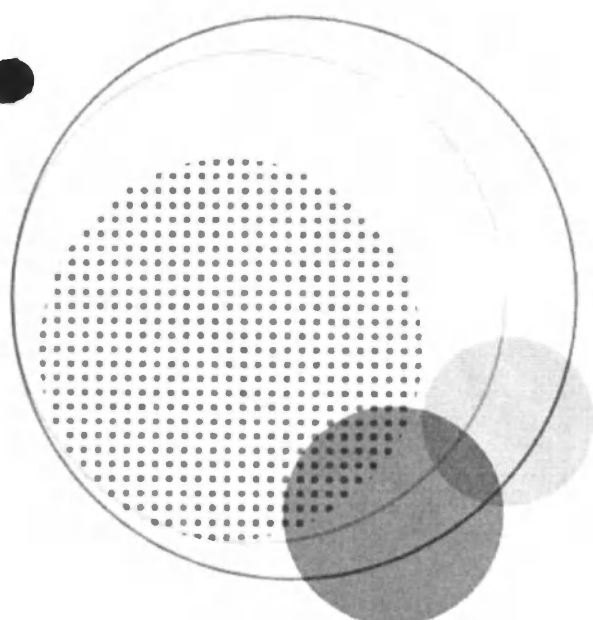
- The use of electric equipment involves electrical risk. The staff has to be capable and have the respective information about the inherent risks of electrical equipment.
- An abnormally leakage of current is symptom of a AC power cable and of power supply unit. If the leakage is over 500 microamperes, personal injuries can be caused.
- Unplug the bed before any bed cleaning or maintenance action is made. Otherwise, personal injuries or material damages may occur. Refer to the specific manual section for further warnings.
- For patient or staff safety, take precaution whilst connecting other wired devices. Avoid harming or pressing the wires of other equipment with the different elements and mobile parts of the bed.
- Concerning the bed connection, avoid harming or pulling the power cord, because it could imply electrical risk.

Entrapment warnings

- There are entrapments or crushing risks whilst lowering the moving parts of the bed base. Check there are not obstacles on its way down.
- There are entrapments or crushing risk between the rail and other bed parts. Check those whilst going down there are no obstacles on its way.
- There are entrapment or crushing risk between the bed parts and accessories (for example the infusion holder) if they are not used according to the manual instructions.

Potential consumers

Trained medical personnel and, to a limited extent, ordinary citizens on the hospital grounds.

Regulations, Rules and Codes

Regulations, rules and codes

NewCare bed is in accordance with the following regulations:

- Regulation (UE) 2017/745
- Directive 2006/42/CE
- Directive 2006/95/CE
- Directive 2004/108/CE

NewCare bed complies with the following rules:

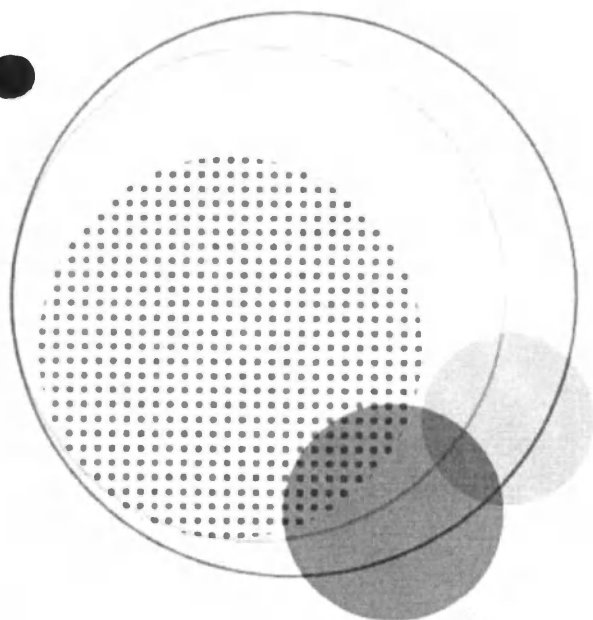
- Standard UNE-EN ISO 12100
- Standard UNE 20324
- Standard UNE-EN 60204-1
- Standard UNE-EN 60601-1
- Standard UNE-EN 60601-1-2
- Standard EN 60601-2-52

NewCare bed has the following marking:

Industrias Pardo s.l. company has the following certifications:

Description

Description

Bed description

The hospitalization bed NewCare is specially created to be used in hospitals, clinics and other similar uses. Designed for 10 to 15 years lifetime, according to its use.

Its structure is made of steel with textured coating painting with high quality polyester base. Thickness of the steel layers between 1.5 and 3 mm (according to zones). This coating is resistant to chemical and mechanical agents.

The bed base is divided in three articulated panels, two of them electrically actioned (backrest and extremities) and another manually articulated (feet). The central panel is fix. The bed base is made of three HPL (High Pressure Laminated) high resistance pieces, removable and fireproof. The one fix panel is made of steel with textured coating painting with polyester.

It is designed to place the patient in different positions (horizontal, fowler, Trendelenburg, etc) with other available emergency pre-programmed positions.

It has four supports for infusion stand, trapeze, Balkan frame, etc,... made of thermoplastic material, as well as drainage and nasogastric supports at both bed sides.

It has four Ø 150 mm wheels. The 5th wheel acts as directional. In case there has not 5th wheel, one of the four wheels of the bed will be antistatic and directional.

Multifunction pedal with three positions:

- Directional wheel activation.
- Wheel releasing and beds free moving.
- Wheels blocking.

Compass lifting and lowering system with great stability.

The backrest and legs modules synchronized motion, whilst activated, improves the pressures distribution and contributes to the reduction of the PU's formation. It has bed base auto regression.

Universal feet module positioning by gas spring.

Trendelenburg and Reverse Trendelenburg positioning system with tilting possibility in both senses.

Different bed movement controls are available, such as, wired patient control, membrane control fitted in the side rails and allows the patient and the caregiver control, and the nurse control that has the option to block the beds movements and has the pre-programmed specific functions.

The bed can be used in 1, 2, 3 y 5 environments according to directive EN 60601-2-52:2010.

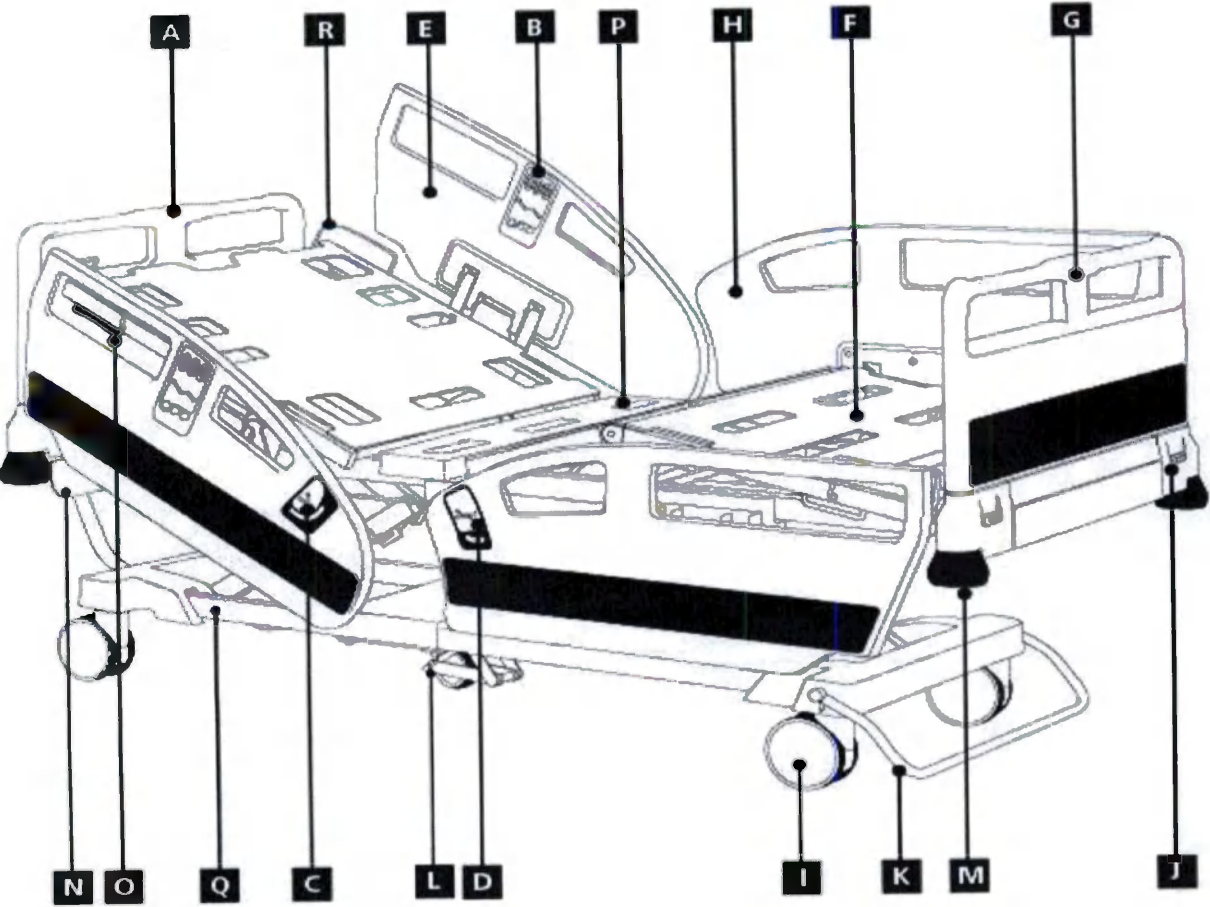
Before using the bed ensure to read and understand the NewCare user manual content. It is important to read and strictly follow the safety instructions contained in the manual.

Intended use

Functional electric 4-section NewCare bed in versions with accessories is designed to care for the patient, increase comfort and create the desired positioning of the patient in a hospital environment .

Potential consumers

Trained medical personnel and, to a limited extent, ordinary citizens on the hospital grounds.

Item	Description	Item	Description
A	Removable curved headboard	I	Double wheels (Ø150 mm.)
B	Membrane control in side rails		Headboard and Footboard unlocking
C	Backrest tilting indicator	K	Braking pedal
D	Trend / Reverse Trend tilting indicator		5 th wheel
E	Backrest side rails	M	Protection discs
	HPL mattress platform	N	Anchored headboard
G	Removable curved footboard	O	Emergency CPR
H	Legrest side rails	P	Sheet steel central section
Q	Bed frame	R	Mattress holders

Common bed weight– 125 kg

Bed weight with accessories -150 kg

Rated power consumption – max 3.15 A (max 400 W)

Safe working load of bed extension– 140 kg

Safe working load drainage holders- 2 kg

Bed components

Electrical components

Do not hit the controls. Do not manipulate in case of breakage and do not access to the controls inside.
In case of power failure, the control will enter in lock mode. To unlock follow the instructions specified in this manual.
In case of breakdown, please, consult your dealer.

The use of electrical equipment carries electric risks. The staff should be capable and must have available the information inherent to the electrical equipment.
The electric circuit is electrically insulated from the bed metal parts.
To avoid motor overload, the manoeuvre should not last over two minutes.
Any electric component reparations should be done with the bed unplugged of mains and disconnected battery.
Do not manipulate the boxes, controls or electrical motors in case of breakdown or access inside the electric components.
Observe that the power socket and bed plug are properly connected. An error respecting this can cause overheating in these components.
Check that the mains voltage match with the nominal voltage needed.
The electrical installation where the bed is connected must have the regulatory protections.

Controls

NewCare bed is equipped with different controls depending on options and configurations that the customer needing. All the controls share the following features:

- Include locking position indicator.
- Compatible with right/left-handed.
- The controls with wired include a hook on rear zone to hang up.

Control box

The control box of the electric system is located next to the extremities motor.

Features:

- Supply voltage: 230 V. 50Hz.
- Operating temperature range: 5° - 40°
- IPX6 protection

Power cord

The power cable has a minimum length of 2.5 m from the mains plug. The power cord has a molded plug. Insulated and sheathed power cables and other external bed flex cables.

Battery (optional)

The battery is located in the bed chassis, in the central area next to the extremities motor and control box.

- Rechargeable battery type: Lead gel.
- Voltage: 24 V
- Capacity: 1,2 Ah.
- Cycle: Approximately 300 movements.

The maximum time a bed with battery can be unplugged is three months, although is recommendable charge it every two months if the bed goes to remain unplugged.

The batteries have a limited life so that it's convenient checking usually its charge condition. It's recommendable replacing them when exhaustion is observed or every three four years.

The nurse control has a battery indicator, which shows when the bed is plugged to mains and the battery is properly charging.

Removing the battery: Unscrew bracket ends and removing bolts which are fastening it. Placing the battery: screw on the bracket ends with bed frame by means of nuts and bolts.

Linear electric actuators

NewCare bed has 2 linear actuators located in bed chassis which take over with the backrest and extremities sections rising; and 2 actuators in the chart for its lifting and lowering.

- Motor voltage: 24 V
- Installation dimensions:
 - 445 mm. (Backrest).
 - 336 mm. (Legrest).
 - 363 mm. (Bed).
- Thrust:
 - 4000 N. (Backrest).
 - 3000 N. (Legrest).
 - 6000 N. (Bed).
- Stroke length:
 - 200 mm. (Backrest).
 - 100 mm. (Legrest).
 - 150 mm. (Bed).
- Protection: IPX6.

The backrest motor includes a clutch for the lowering in case of emergency and self-locking to avoid entrapment.

Mechanical components

Side rails

The use of any element to cover the mattress reduces the rails real height as well as other equipment dimensions. In case of using such elements, evaluate the patient falling or entrapment risks and take the necessary measures.

When the bed is moved, with or without patient in its lowest position, the rails must be in its highest position to avoid patient injuries whilst moving.

Installed in the bed sides, the objective of the rails is to remember the patient the bed limits but it will never be a contention system it selves.

Plastic side rails

There are two folding rails at the backrest section and articulates with it (with optional control membranes in the inside and outside of the rail) and two folding rails fixed to the bed chassis at the legs section.

Rails made of high resistance polyolefins, installed on sides of the bed and fixed to the bed base of bed.

Each side rail is provided with a lever to be folded with one-hand operating system. Its drop is cushioned with a gas spring which avoids sudden fall. Side rails have a safety clamp in highest position and anti-trapping system which avoids they been locked in lowest position.

Each rail has a lever for folding and a degrees indicator that, in case of the legs section rails. Indicates the degrees of trendelenburg and reverse trendelenburg positions and in case of the backrest rails shows its angle of inclination.

Metal rails in stainless or painted (optional):

Installed on sides of bed. 4 tubes side rails in stainless-steel or powder coated assembled on close and rounded tube brackets made in injected polyamide with steel core.

Each side rail is provided with locking/unlocking system.

Bed extension

When the bed is on its extended position, do not seat or lean on the feet side because it might overturned the bed.

The bed is 300 mm extensible, to accommodate the bed to the physical patient characteristics.

Pneumatic adjustment of the leg section of the bed

Raise the leg section only when the patient is in a bent knee position, otherwise there is a risk of injury

The leg section of the stock can be raised or lowered to the desired position using a pneumatic spring. Changing the position of the leg section of the bed helps patients with medium and long stay to change position, and also helps prevent the formation of pressure ulcers in the ankle area.

Feet section position

Articulate only the feet module when the knees are bent, otherwise it could cause patient injuries.

The feet module can be raised or lowered by a gas spring until the chosen position is reached.

This positions helps the medium and long stay patients to change positions and contribute PU's prevention in the heels areas.

Headboards and footboards

It is a cart fixed headboard system that maintains it always in the same position and height, independent to the bed positions and moving made.

Removable head and foot board with locking system

Located in the head and feet sides of the bed, its objective is to remember the patient the bed limit but they are not a contention device it selves.

Removable foot/headboard in polypropylene. They are designed with rounded edges and ergonomic grips which improves handling during transport.

With a bed lock/unlock system in the centre of them that makes it possible to place and remove easily by the bed operator without any tools.

Wheels

NewCare beds has **four integral wheels** (Ø 150 mm.), with high rolling, wear and rust resistance and it has a minimum noise emission whilst moving.

Optionally, it is possible installed different wheel models: Double or metallic wheels both, with 150 mm of diameter.

Braking of bed

Remember always to keep the brake on when the bed is not being moved, to avoid injuries in case of accidental displacements.

Ensure, at any bed transportation, that it is unplugged from mains

Multifunction bar for centralised braking

Placed at the feet side of the bed, the multifunction bar allows the centralised braking with auto locking of the four wheels, total **free movement** and a third possibility, depending on the bed configuration the 5th wheel can be activated or the **directional or 5th wheel** instead, once they are activated the beds manoeuvrability will be improved whilst it is being moved, avoiding lateral displacements.

Lateral pedals for centralised braking (optional)

Placed next to the wheels of footboard or the four wheels. These pedals allow the centralised braking with auto locking on the four wheels, the total **free of movements** and the **directional or 5th wheel** activation.

5th wheel (optional):

It has an unrestrained Acoustic Warning as an option.

Made of high resistance synthetic plastic material to avoid exposed metal parts and, therefore, making impossible the danger of corrosion.

Some of its characteristics are:

- **Cushioned:** The cushioning system avoids the danger against the presence of small steps.
- **Antistatic:** It allows the electricity conduction to the ground because of the wheel material composition.
- **Directional:** It locks the fifth wheel in directional position, relieving the sanitary staff of extra effort because it facilitates the manoeuvrability of the bed.

RCP bilateral actioning system

Do not use the RCP bilateral actioning levers to lift the beds backrest under no circumstances, because it may cause a lifting motor failure and the bed will stop working.

Placed in both sides of the back of bakrest section, the CPR bilateral actioning system allows the lowering of the bed quickly.

Accessories support

These supports can only be used to hold our manufacturing products, otherwise, product damages or patient injuries may be caused.

Along all the bed different supports can be found to fit accessories, such as:

- Corner supports to fit infusion stands, trapezes or Balkan frames.
- Drainage holder arcs under the bed base, with regulated hooks.
- Patient fixing arcs with “SEGUFIX” system, placed along the sections.

Protection discs

Located in the four bed corners, the protection discs protect the bed against shocks.

NewCare bed disposes of two models of discs, according the external dimensions that the customer required.

Bed model	Color	diameter	picture
Standard	Blue (wide)	Ø104 mm	
S	White (narrow)	Ø84 mm	

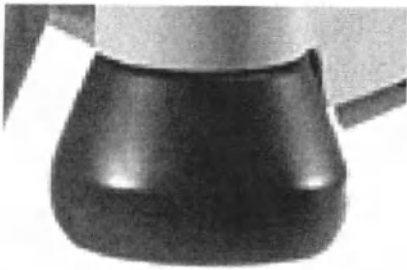
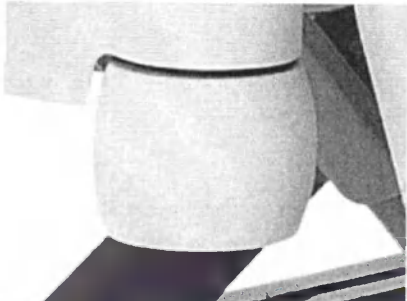
XS	Blue (wide)	Ø104 mm	
XSS	White (narrow)	Ø84 mm	

Table under footboard

Under the bed base at the feet section, the bed has a high resistance extendible tray that can be used by the patient or the assistance personnel.

Under it you can find the nurse advanced control, in case it is included.

SWL 5 kg

Multifunction table

Multifunction table placed in footboard using non-tooling hooks which provide good stability and reliability. Could be used by healthcare personnel or by patients when eating or other activities.

SWL 5 kg

Rx cassette tray

Thanks to guides under backrest, it's possible including a X-Ray cassette holder. This tray is made of stainless-steel with upper handle which supports during traction, and it's high adjustable to adapt any patient. X-Ray cassette tray is removable from upper section without removing mattress base

Light under bed

Under frame light to aid healthcare personnel during transport through dark corridors or patient for bed exit during the night.

The light can be turned on anytime in side rails membranes with this option as well as in the patient control with light option.

Light turns on automatically when the bed is in the lowest position.

Brake alarm

NewCare bed can be equipped with brake alarm. This acoustic alarm is activated when the bed is plugged and the brake released.

Accessories

Check regularly the correct state of accessories.

IV-Pole

The infusion stands are located in one of the supports and it is made of two painted steel or stainless steel parts, a handle to regulation height and fixing label with liquid protection system in the lower tube.

Two or ~~four~~ holding hooks system with 2 kilograms capacity per hook.

Trapeze bar non-adjusting

The trapeze bar is located in one of the supports and it is composed of one steel arm with ergonomic designed handle and a belt with anti-sliding system. 300 mm length belt, non adjusting

Trapeze bar adjusting

The trapeze bar is located in one of the supports and it is composed of one steel arm with ergonomic designed handle and a belt with anti-sliding system. 300 mm height adjusting system, by pressing button.

Balkan frame

Aluminium frame structure with octagonal section.

It is composed by two curved arcs fixed on both sides of headboard and footboard and one central bar along the length of bed

The kit is composed by few jaws, balancers, trapeze (with handle and belt) and infusion stand.

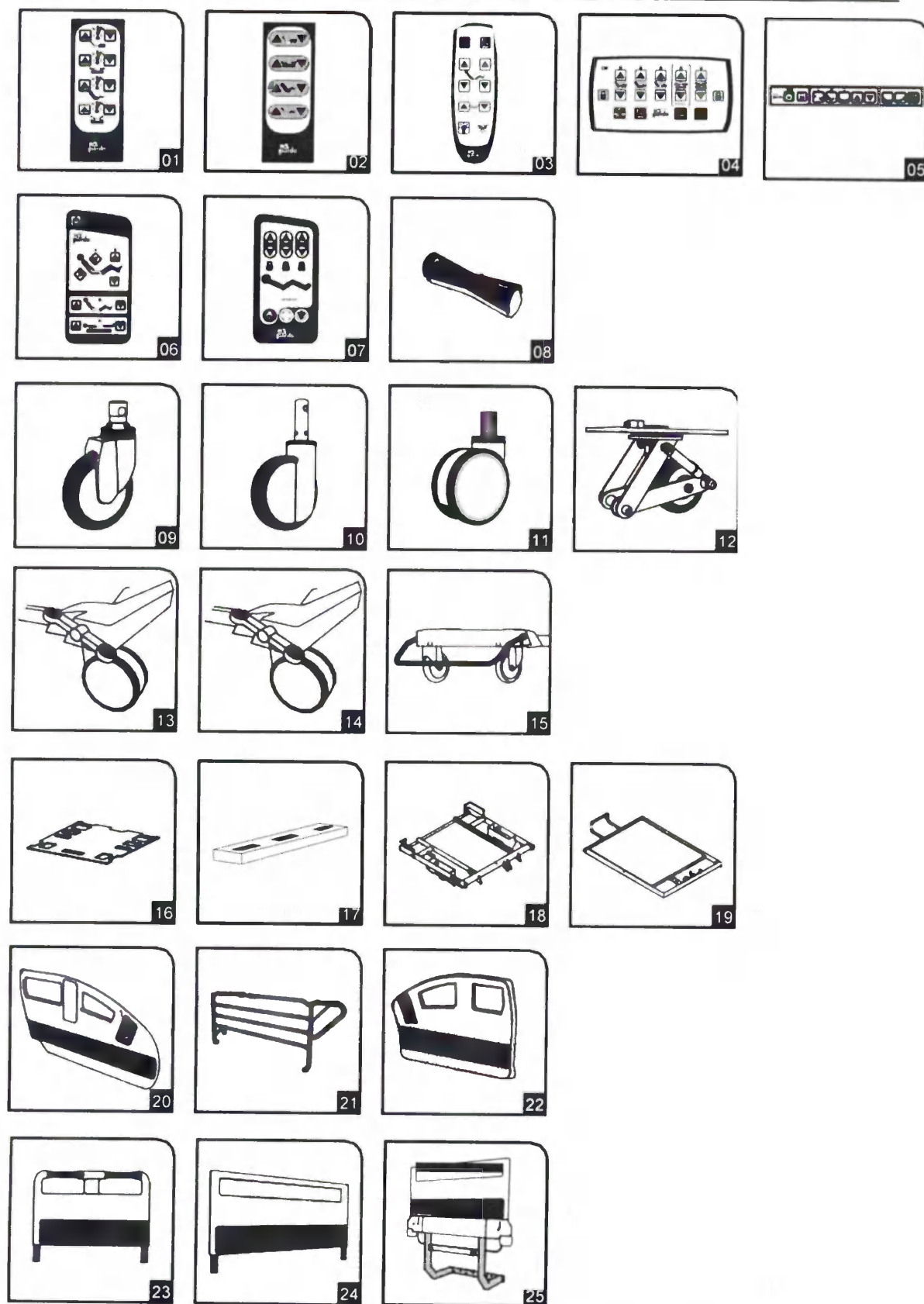
SWL 75 kg

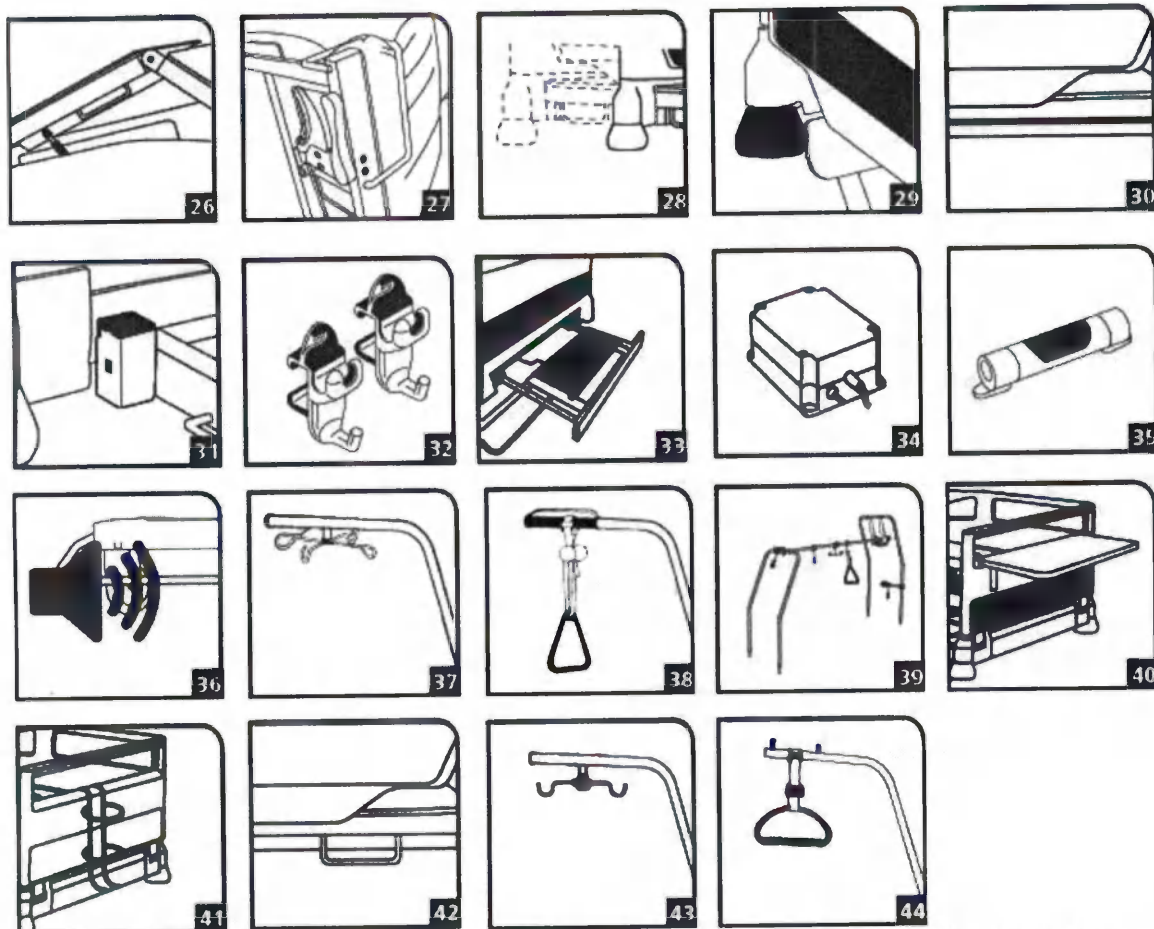
Tank holder

Metallic structure where can be attached to the footboard or the headboard of the bed and can hold an oxygen bottle.

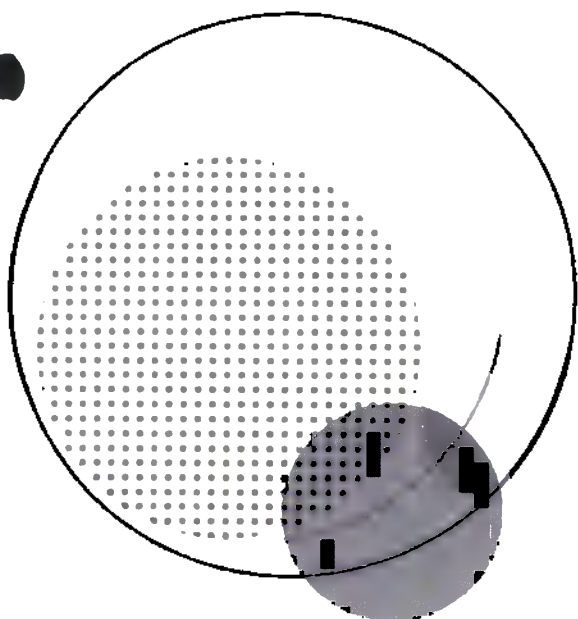
SWL 15 kg

List of components

- | | | |
|---|---------------------------------------|------------------------------|
| 01 Patient hand control with locking c | 16 HPL mattress platform | 31 Accessories supports |
| 02 Patient hand control without locking | 17 Sheet steel central section | 32 Drainage holder |
| 03 Patient hand control with light | 18 Guides for X-Ray tray | 33 Table under footboard |
| 04 Nursing hand control | 19 X-Ray cassette tray | 34 Battery |
| 05 Nursing control on footboard | 20 Integral side rails | 35 Light under bed |
| 06 Membrane of side rails | 21 Four tube side rails | 36 Braking alarm |
| 07 Membrane with light of side rails | 22 Legrest rail for safe exit | 37 IV-Pole |
| 08 Lifting pedal | 23 Removable curved head/foot board | 38 Trapeze bar non adjusting |
| 09 Metallic wheel Ø150 mm. | 24 Removable straight head/foot board | 39 Balkan frame |
| 10 Integral wheel Ø150mm. | 25 Anchored headboard on cart | 40 Multifunction table |
| 11 Double wheel Ø150mm. | 26 \Spring gas for feet | 41 Tank holder |
| 12 5 th wheel | 27 Emergency CPR | 42 Patient strap holder |
| 13 2 centralised braking pedals | 28 Extensible bed base | 43 II pole |
| 14 4 centralised braking pedals | 29 Protection disc | 44 Trapeze bar adjusting |
| 15 Centralised braking bar | 30 Fixing arcs | |

Bed operation

Scope of delivery

Reception

- When the product is received, check if the shipment is correct.
- Inform the forwarder and supplier when any deficiency or damages are detected immediately.
- Unpack the product conveniently and eliminate the packaging waste as it is explained further down.

Installation

- Ensure not to step on any wire whilst moving the bed.
- Ensure the power supply cable is attached to the headboard hook.
- Ensure the wheels are unlocked before moving the bed.
- Move the bed on appropriate surfaces.
- Mount the headboard and footboard on the bed.
- Ensure the power supply cable is not trapped or overstretched when the bed is placed.
- Ensure the plug is correctly connected

Initial operation

- Our beds are prepared to operate from the outset. The battery is charged, in case, at the timethe bed is plugged, the battery is discharged, it must be recharged before putting it into service.
- An initial bed movement testing can be made:
 - Lift the bed to its highest position.
 - Take down the bed to its lowest position.
 - Articulate the bed base sections.
- Check the rails cushioned lowering (according the configuration of bed)
- Check that the headboard and footboard can be released correctly.
- Check the correct function of all the bed controls

*Only authorized maintenance staff of the establishment should do the star up of the beds.
The manufacturer will not be responsible of any bed or battery damage caused because by an inappropriate handling of the bed or components*

Controls functioning

The bed runs by pressing different control elements. The control elements will depend on the chosen bed model.

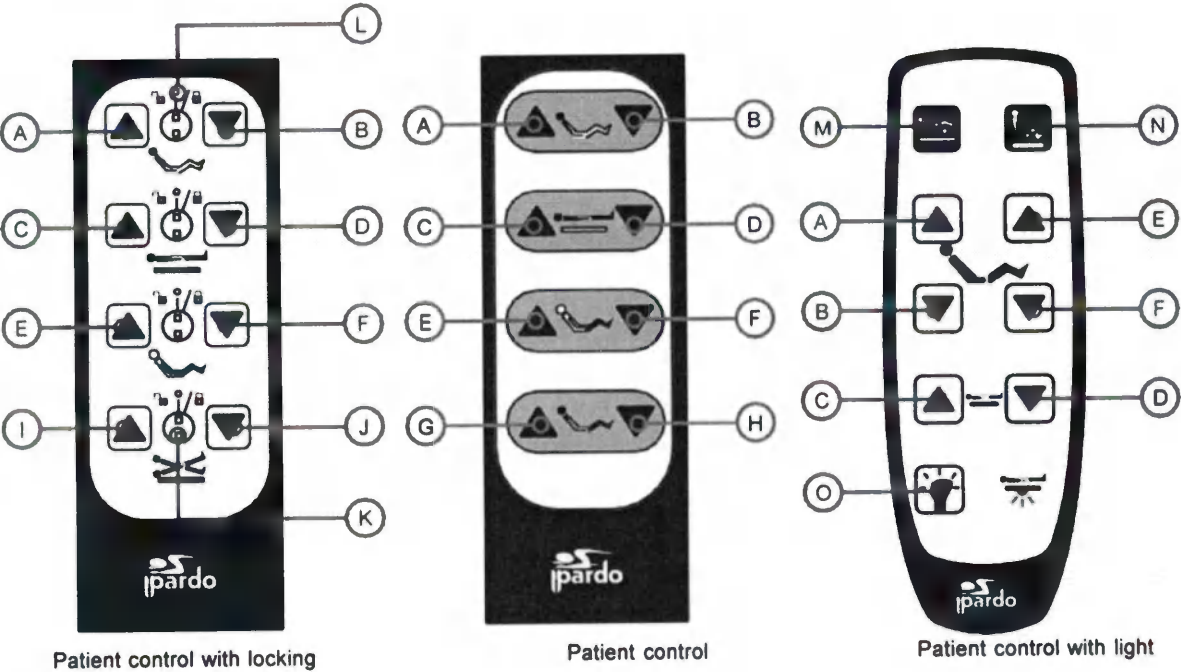
The individual function locking in nurse control affects to all the other control elements.

Kind of controls

Wired hand control:

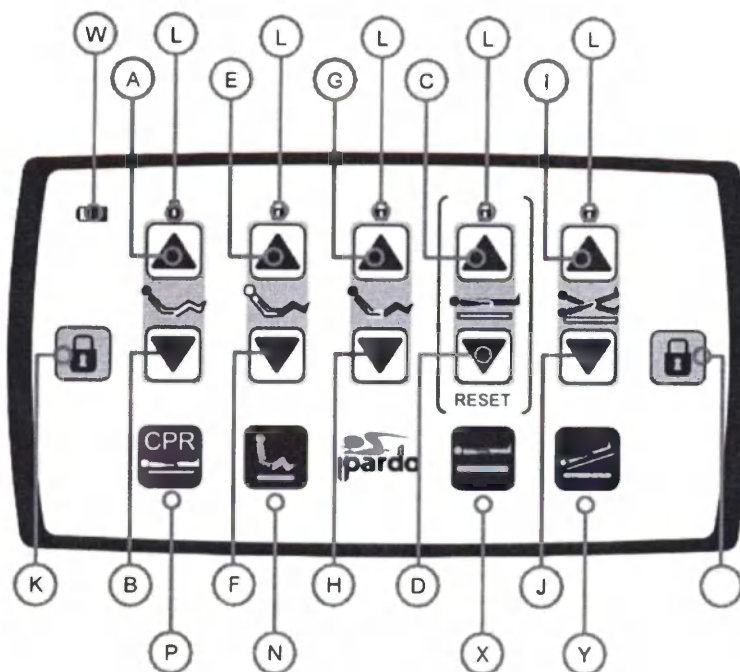
NewCare bed has three models:

- Patient control with locking.
- Patient control without locking.
- Patient control with light.

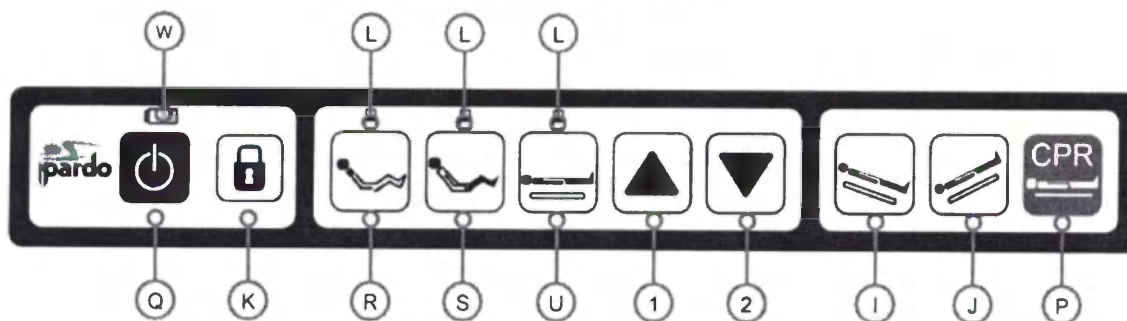
A	Rise backrest	I	
B	Lower backrest	J	Reverse trendelenburg
C		K	Locking movements
D	Lower mattress platform	L	Indicator Led of locked movements
E		M	Flat position (Lower all sections simultaneously)
F	Lower legrest	N	Cardiac chair position
G	Rise backrest and legrest simultaneously (autocontour)	O	Light under
H	Lower backrest and legrest simultaneously (autocontour)		

Nursing controls:

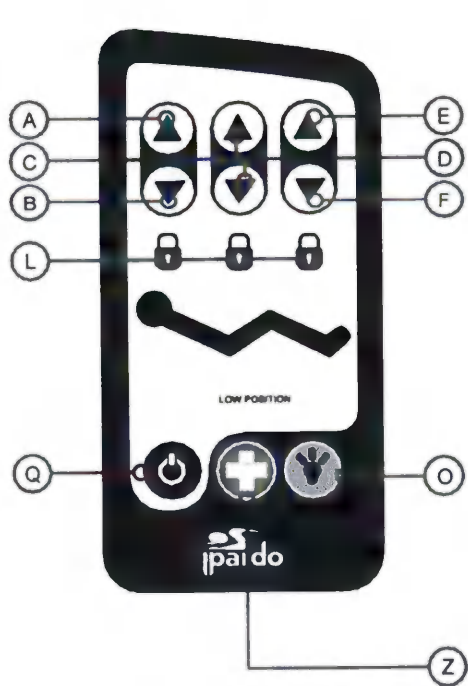
Nursing control

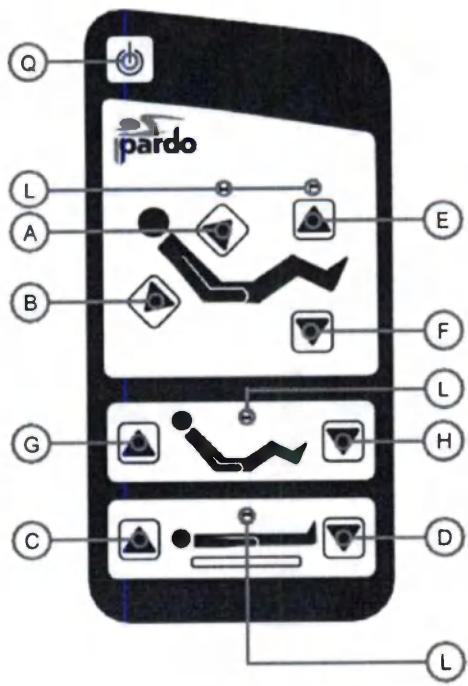
Membrane on footboard

A	Indicator Led of locked movements	N	Cardiac chair position
B	Lower backrest	O	Light under bed
C	Rise mattress platform	P	Emergency CPR
D	Lower mattress platform	Q	On / Off button (to active the membranes)
E	Indicator Led of locked movements	R	Backrest selector
F	Lower legrest	S	Legrest selector
G	Rise backrest and legrest simultaneously (<i>autocontour</i>)	U	Bed base selector
H	Lower backrest and legrest simultaneously (<i>autocontour</i>)	W	Battery indicator
I	Trendelenburg	Y	Emergency trendelenburg position
J	Reverse Trendelenburg	X	Explored position (maximum height)
K	Locking movements	1	Higher
L	Indicator Led of locked movements	2	Lower

Membrane controls:

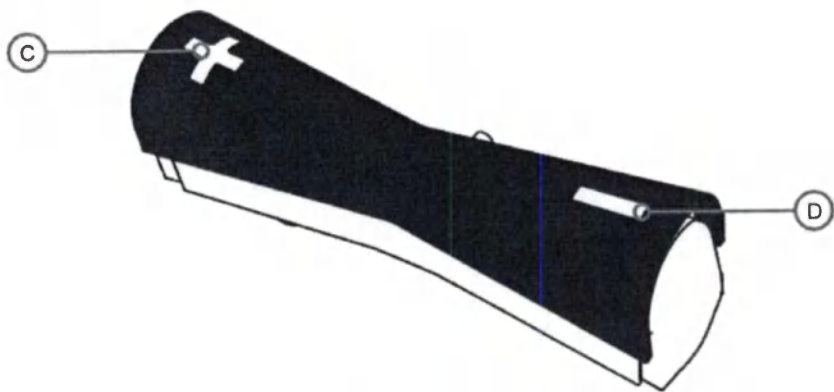
Membrane control with light

Membrane control

A	Rise backrest	G	Rise backrest and legrest simultaneously (<i>autocontour</i>)
B	Lower backrest	H	Lower backrest and legrest simultaneously (<i>autocontour</i>)
C	Rise mattress platform	L	Indicator Led of locked movements
D	Lower mattress platform	O	Light under bed
E	Rise legrest	Q	On / Off button
F	Lower legrest	Z	Nursing call

Lifting pedal:

C	Rise mattress platform	D	Lower mattress platform
---	------------------------	---	-------------------------

Operating

Basic control operating

Do not seat on the articulated sections whilst raised, neither activate the movements if there are people seated in those modules.

There are entrapments or crushing risks whilst lowering the moving parts of the bed base. Check there are not obstacles on its way down

There are crushing risks whilst lowering the bed between it and the cart, in case you are resting your feet on the cart.
There are crushing risks whilst lowering the bed when the rails are all the way down, or if when the rails are being folded with the bed in its lowest position.

Rise and low backrest

With the backrest rising and lowering buttons (A and B) the backrest module can be positioned at the desired angle

To activate the movement keep on pressing the backrest rising or lowering button to reach the maximum or minimum tilting respectively.

From footboard membrane it's possible doing this action selecting backrest button (R) and pressing raising or lowering arrows (1 and 2).

Do not use the RCP bilateral actioning levers to lift the beds backrest under no circumstances, because it may cause a lifting motor failure and the bed will stop working.

Rise and low legrest

With the backrest rising and lowering buttons (E and F) the legrest module can be positioned at the desired angle

To activate the movement keep on pressing the legs raising or lowering button to reach the maximum or minimum tilting respectively.

From footboard membrane it's possible doing this action selecting legs button (S) and pressing raising or lowering arrows (1 and 2).

Rise and low mattress platform

With the bed rising and lowering buttons (C and D) the bed can be positioned at the desired height.

To activate the movement keep on pressing the bed rising or lowering button to reach the maximum or minimum height respectively.

From footboard membrane it's possible doing this action selecting bed height button (U) and pressing raising or lowering arrows (1 and 2).

The maximum height improves the assistant ergonomics and the minimum height the patients.

Do not try to lift or low the bed with other means than the ones designed by the manufacturer. If this warning is not taken in count a lifting mechanism dysfunction or body injuries may occur.

*There are crushing risks whilst lowering the bed between it and the cart, in case you are resting your feet on the cart.
There are crushing risks whilst lowering the bed when the rails are all the way down, or if when the rails are being folded with the bed in its lowest position.*

Rise or low the backrest and legs simultaneously (autocontour)

With the backrest and legs rising and lowering buttons (G and H) the backrest and legs modules can be positioned at the desired angle at the same time.

To activate the movement keep on pressing the backrest and legs rising or lowering button to reach the maximum or minimum tilting respectively

Trendelenburg & Reverse trendelenburg

With the Trendelenburg (J) and reverse Trendelenburg (I) button:

- Trendelenburg

Activation: Press the Trend (I) button to tilt the footboard related to the headboard until the maximum inclination is reached. the rest of the bed base stays in the original position.

- Reverse trendelenburg

Activation: Press the reverse Trend (J) button to tilt the headboard related to the footboard until the maximum inclination is reached

Trendelenburg and reverse Trendelenburg movements always have to be done under medical supervision.

Pre-programmed functions

Emergency CPR

Through CPR button (P), all the bed modules rest in horizontal position and the bed will go down to its lowest position, for a possible urgent cardiopulmonary resuscitation.

There are entrapments or crushing risks whilst lowering the moving parts of the bed base. Check there are not obstacles on its way down.

Explored position (maximum height)

Via maximum height button (X) the bed reaches its highest and horizontal position, to facilitate a doctor's ergonomic position

There are entrapments or crushing risks whilst lowering the moving parts of the bed base. Check there are not obstacles on its way down.

Emergency Trendelenburg position

Activation: Press the emergency Trendelenburg button (Y) to place the bed base sections horizontal and tilting the footboard towards the headboard.

Trendelenburg movement always has to be done under medical supervision.

Flat

Position that can be selected only from the patient control with which, by pressing the pertinent button (M) the bed's body sections go down to their lowest position leaving the bed in a totally horizontal position what helps to the comfortability and wellness of the patient whilst resting.

Cardiac chair position

With the cardiac chair button (N) the modules tilt to its maximum; the bed base tilts to the reverse trendelenburg position. The movement successions needed to reach this position follows certain order and intervals to improve patient ergonomics and safety.

There are entrapments or crushing risks whilst lowering the moving parts of the bed base. Check there are not obstacles on its way down.

On / Off button

By pressing the On / Off (Q) button, the rails membrane control is activated to be able to use the functions located in it. After an inactivity of 60 seconds, the membrane turns off and it is necessary to press again the button to reactivate.

Each membrane has its own activating button. It is necessary to press its own button to activate each membrane.

Nursing call button (needs a H.I.S. installation)

Nurse call system needs H.I.S installation (Hospital Information System) in the hospital which consists in software managing each variable in hospital facilities, personnel and patients.

In case of H.I.S. is available, membranes embedded in side rails will be able to connect to the system.

Pressing (Z) button nursing call will be activated in case of patient or relatives need assistance. To active it just press the button once.

Light under bed button

The under bed lighting button (O) lights on the luminary located under the bed base and allows to walk by night without turning on the rooms lighting, without bothering near people.

To turn on the light press the button once and to turn it off press it once again.

Light indicators

Movements lock indicator

If the indicator (L) is lighted, shows that the corresponding movement is locked.

Low battery indicator

The nursing controls have a light indicator (W) to remark if the battery is charging. When the LED is lighted the battery is charging, if it is off, the battery is completely charged.

Low Position indicator (only in membranes with light)

The membranes with light have a low position indicator LED.

The LED is automatically lighted when the bed is on minimum height position.

Bed function locking

Patient control with locking

Function locking is used to disable bed basic functions, mainly to lock movements which has to be done under medical guidance.

Locking: introduce locking key in the central zone of the function to be disabled and turn it clockwise until yellow indicator is displayed (blue circle in attached picture).

Unlocking: turn the key counter-clockwise until green indicator is displayed.

Nursing control

Locking different sections or functions will be made pressing the button signalled with a key (K) and then press at the same time the section to be locked (A, C, E, G or I). Indicator light (L) which has been pressed will become fixed. To unlock, redo the same action.

The indicator light will remain flickering for a few seconds.

Nursing membranes on footboard

- Backrest movement locking

Press the indicated locking button with a key (K) and then at the same time press backrest button (R). The indicator light (L) will remain fixed. To unlock, redo the same action. The

indicator light will remain flickering for a few seconds.

- Legs movement locking

Press the indicated locking button with a key (K) and then at the same time press legrest button (S). The indicator light (L) will remain fixed. To unlock, redo the same action. The

indicator light will remain flickering for a few seconds.

- Raising/lowering movement locking

Press the indicated locking button with a key (K) and then at the same time press platform mattress button (U). The indicator light (L) will remain fixed. To unlock, redo the same action. The

indicator light will remain flickering for a few seconds.

Other control functions

System reset

System reset is necessary when all system is locked due to a security limits has been exceeded. A beep will sound from the control box warning there has been a system failure.

- Patient control

Keep pressed for some seconds buttons C & D until the beep disappears. The system will reboot and it will work properly again.

- Nursing control

Keep pressed for some seconds buttons C & D (the buttons are marked) to reset the system for a few seconds until the beep disappears. The system will reboot and it will work properly again

- Nursing membrane on footboard

If it's necessary doing system reset from membrane control placed in footboard, this is done pressing buttons 1 & 2 at the same time up to the beep disappears.

Bed components functioning

CPR bilateral actioning system

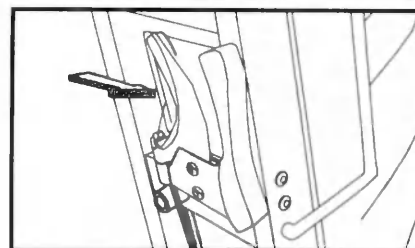
Do not use the CPR bilateral actioning levers to lift the beds backrest under no circumstances, because it may cause a lifting motor failure and the bed will stop working.

There are entrapments or crushing risks whilst lowering the moving parts of the bed base. Check there are not obstacles on its way down.

Activate by hand any of the two levers placed in the beds centre section. For patient safety, the system has a shock absorber to avoid a sharp fall of the backrest.

Rails

When the rail has been raised, ensure it is fully locked because it could produce damages if it accidentally falls down.

Whilst the bed is being moved, with or without patient, in its lowest position, the rails should always been raised because equipment damages or patient injuries may occur.

Do not seat or lean on the rail because it might harm it.

The use of any element to cover the mattress reduces the rails real height as well as other equipment dimensions. In case of using such elements, evaluate the patient falling or entrapment risks and take the necessary measures.

Designed to be used in our manufactured beds.

Regularly check the tight of the bed fixings.

There are crushing risks whilst lowering the bed when the rails are all the way down, or if when the rails are being folded with the bed in its lowest position.

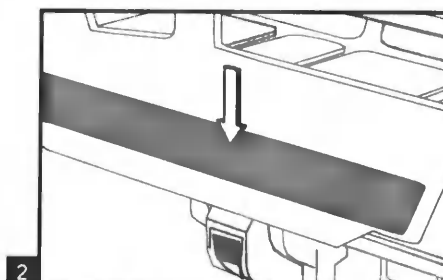
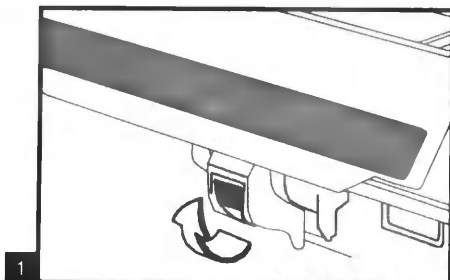
Plastic Full length rails

Rail elevation

To rise a rail you have to pull up until its highest position and lock.

Rail lowering

To lower a rail pull the underneath hanger (1) outwards. The rail (2) will soften come down.

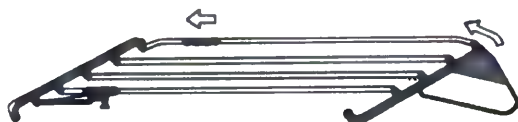
Metal Tube side rails

Side rail raising

Pull of the side rail (1) towards headboard up to hearing “click” from the lower handle (2), then the side rail is locked in highest position.

1

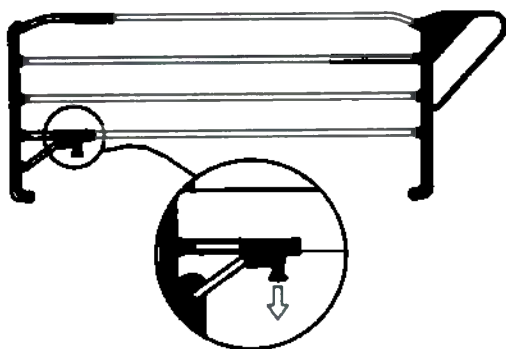
2

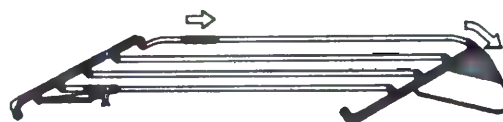
Side rail lowering

- Unlock side rail lower handle pulling of it down (3).
- Fold the said rail towards footboard (4) taking it until lowest position where it remains locked.

3

4

Removable Headboard / Footboard with locking system

Placement

Place the footboard anchors over the beds corner holes.

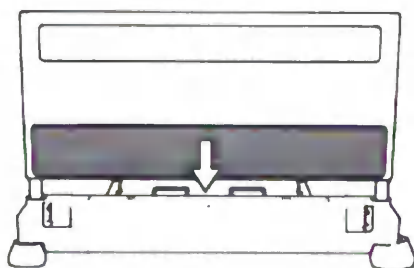
Insert until they are fully in place (1). Lock by turning the levers at each side of the bed towards the inside of the bed until they are horizontal, as shown in the drawing (2).

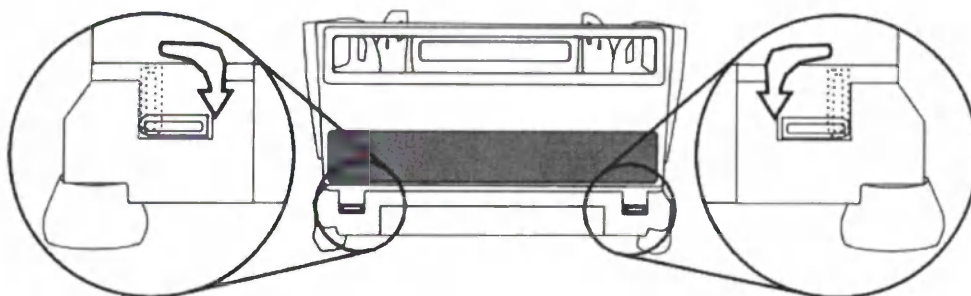
Ensure, by pulling up, that the rails are fully locked.

There are entrapments or crushing risks when the headboard or footboard is being placed. Ensure that there are no obstacles during operation.

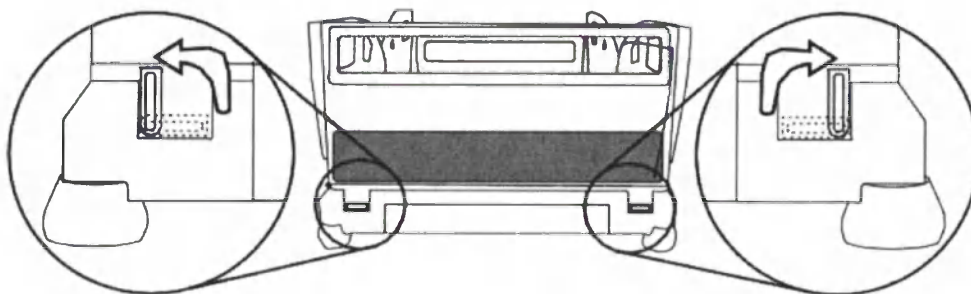
1

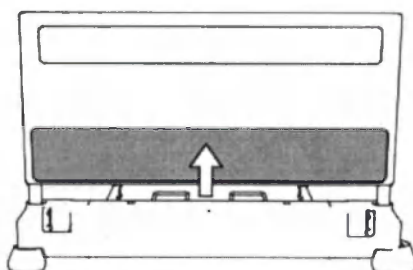
2

Removal

Unlock by turning the levers on both sides (3) of the bed towards the outside of the bed until they are vertical, as shown in the drawing. Hold the foot board and pull it up in a vertical direction (4).

3

4

Multifunction pedal

Remember always to keep the brake on when the bed is not being moved, to avoid injuries in case of accidental displacements.

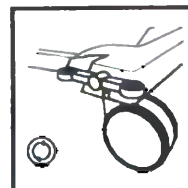
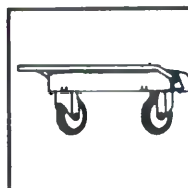
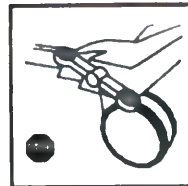
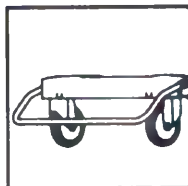
Braking of bed

When the brake is active, the four wheels should be locked simultaneously to keep the bed still.

For actioning step on the pedal all the way down.

This operation is the same for the centralised braking bar and the butterfly pedals located in each wheel.

In case it has a 2 + 2 brake, the pedals only brakes the front or rear wheels.

Moving of bed

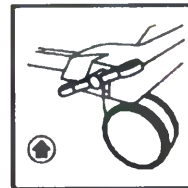
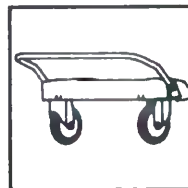
In bed moving position, the wheels remain free to move the bed in any direction.

To activate push the pedal until its intermediate position is reached.

Directional or 5th wheel actioning

The actioned directional wheel allows to carry the bed straight to facilitate the movement without swinging.

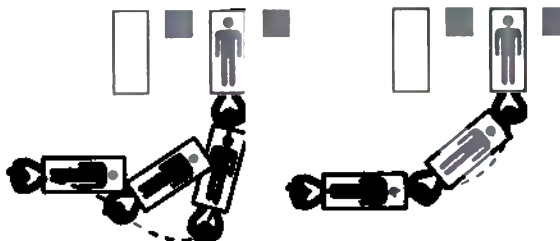
To active it, lift the pedal to its highest position.

*Optionally, the bed has an unlock acoustic, which beeps when the brake pedal is not activated and the bed is still plugged.

5th wheel

The bed can have a 5th wheel as an option. The active 5th wheel makes easier the bed transportation along the corridors (keeping straight). Having a rolling point in the middle of the bed allows to make sharp turns with less effort (entering in an elevator or turning a corner). Another benefit that brings the 5th wheel is to make a complete 360° turn without moving the bed towards any other direction.

Cushioned, antistatic and directional

- The cushion system avoids the damages of the presence of little gaps.
- Allows to conduct the electricity to the ground because of the wheel composition.
- Allows to lock the fifth wheel in directional position, makes it easier the bed manoeuvrability.

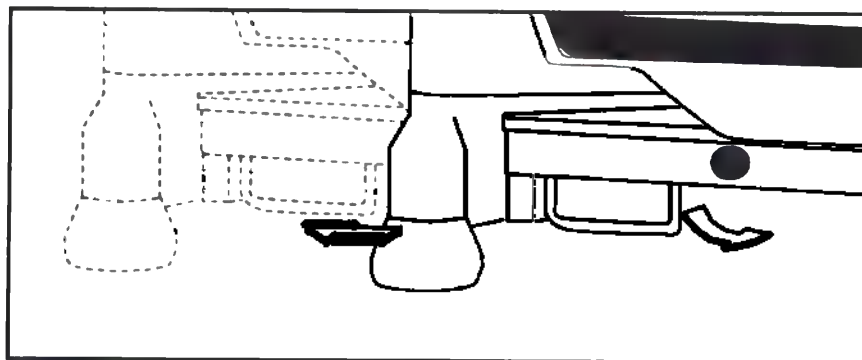
Bed extension

*When the bed is on its extended position, do not seat or lean on the feet side because it might overturned the bed.
Ensure to do the bed extension with the brakes on to avoid involuntary movements.*

Extension of bed

Operate the locking knobs on both sides of the bed by pulling and turning them 180° to keep them unlocked. Take the foot section by the protective discs and pull outwards until it cannot go any further. Lock the knobs on both sides of the bed, turning them 180° again, and pull or push the foot section until it is completely locked on both sides of the bed.

Removing from bed

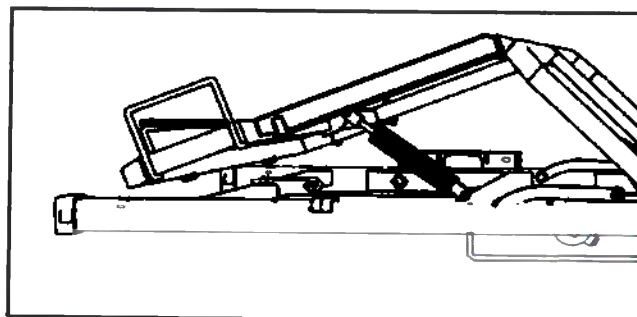
Operate the locking knobs on both sides of the bed by pulling and turning them 180° to keep them unlocked. Take the foot section by the protective discs and push inwards until it cannot go any further. Lock the knobs on both sides of the bed and pull or push the foot section until it is completely locked on both sides of the bed.

Feet module positioning system

It is recommended to articulate the feet module only when the knees are bent, otherwise it could cause patient injuries.

Once the legs module has been electronically raised, move the feet section to reach the desired position. The damper system allows fixing the feet in the desired position. To increase the angle only needs to pull upwards the module. To reduce the tilting, push the module down until it is in the desired position.

Extended table under footboard

Before removing or introducing the table, ensure the brakes are on. Ensure the table is locked in its position when it has been inserted.

Table extraction

To unfold the under footboard table you have to hold it from below and pull outwards until it stops

Table folding

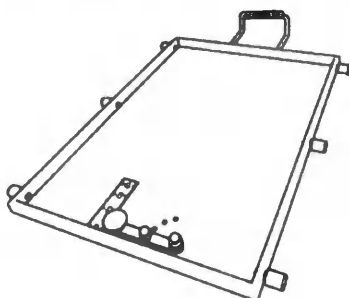
To insert back the table, hold it from below and push inwards the bed until it is locked in its original place.

Cassette carrier

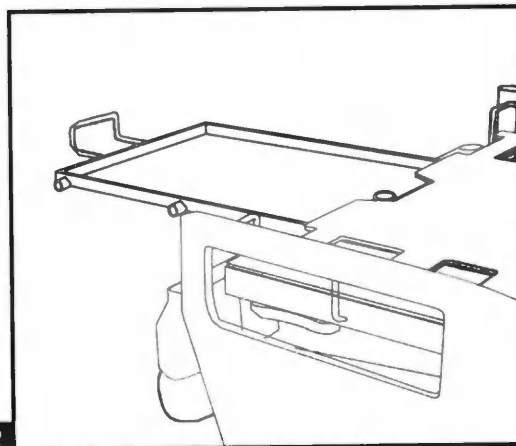
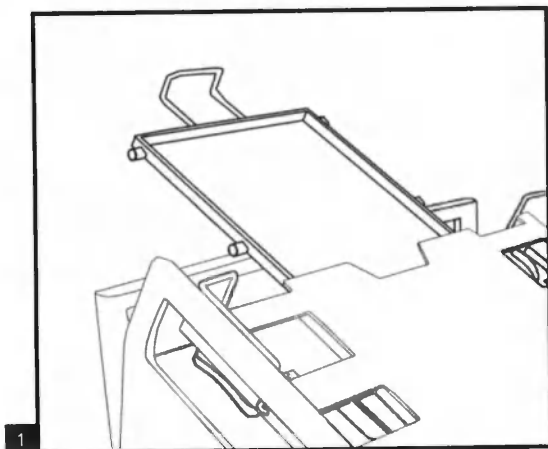
The bed has, as an option, the possibility to include an X-ray cassette carrier to be able to make radiographs to the patient of the thorax section without transferring the patient to a x-ray stretcher or bed.

The cassette carrier has a height adjuster to place the x-ray cassette in the desired place.

To do so, the backrest section has two guides to place the cassette carrier. It can be placed in two ways:

- Rise the backrest until it is over the headboard and be able to introduce it (1).
- Remove the headboard and introduce the cassette carrier without the necessity to lift the backrest avoiding to move the patient (2).

Bed base panels

Ensure by pulling the panels up, that they are correctly installed.

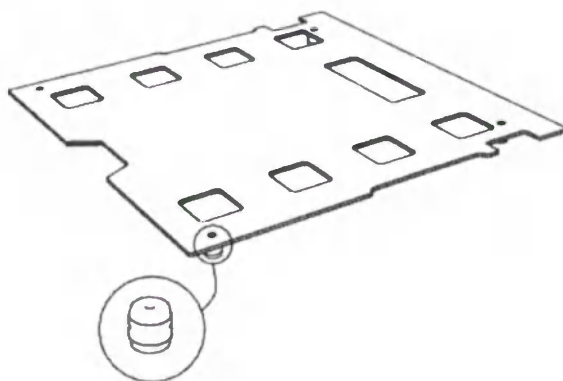
There are entrapments or crushing risks during the bed base panels fixation. Ensure that there are no obstacles during operation

Double regression HPL removable panels, both thorax and legs sections, to prevent PU's.

The thorax section is prepared to include an x-ray cassette carrier as the material the panels are made are radio translucent.

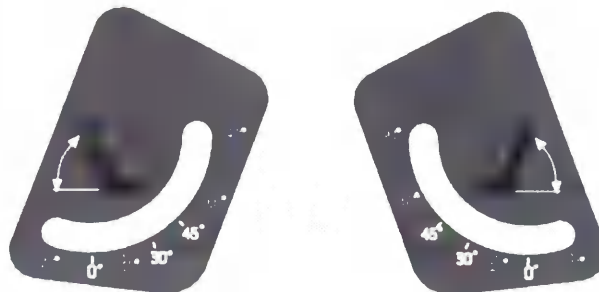
The panels are easy and friendly removable as it does not need to use any tool to fix them. They have a clipping piece to place them on the beds structure.

Movement indicator

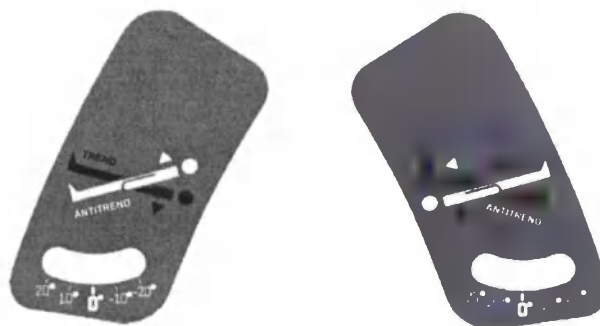
Backrest

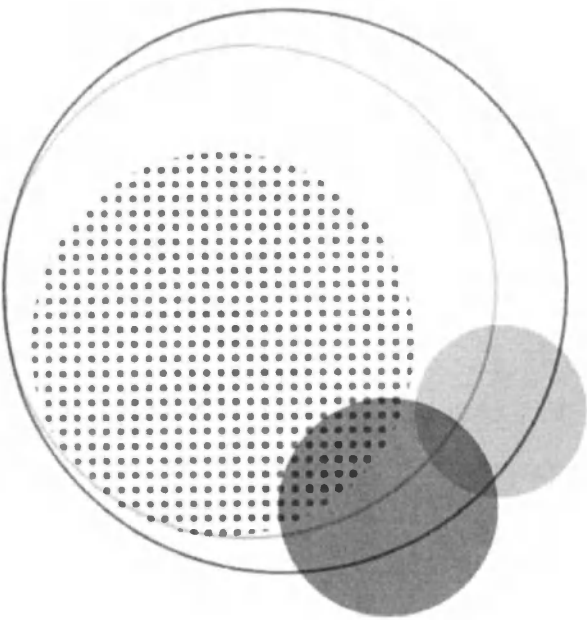
Located in the backrest rails, it shows the tilting degrees in the rising and lowering movements of this section.

Trendelenburg and reverse trendelenburg

Located in the extremities rails, it shows the tilting degrees in the trendelenburg and reverse trendelenburg movements of the bed

Cleaning and maintenance

Bed cleaning

Do not use any steaming cleaning systems because an excessive humidity and temperature may damage some components.
Do not use organic dissolvent for bed cleaning.
Do not use pressure water jet

Any cleaning operation has to be done with the bed unplugged.

Before cleaning:

- Unplug the power cord from the mains
- Move the wheel pedals to the brake position

During cleaning:

- Do not use abrasive cleaners
- Do not use water above 50 ° C
- Do not use steam cleaner, tunnel wash, high pressure hoses
- Do not immerse the equipment in water
- Do not clean or degrease the steel parts of the drives
- Use only approved cleaning agents. If it is necessary to use a disinfectant, strictly follow the instructions for its use.
- Keep sharp objects such as razors, etc. away from the mattress and cover
- Do not use petroleum-based products
- Do not disconnect the control panel or actuators from the control unit when cleaning

GENERAL CLEANING INSTRUCTIONS (ALL SURFACES)

The washing up:

- Wipe with a sponge or soft cloth dampened with warm water containing soap or mild detergents
- The mattress cover is machine washable (95 ° C)

Drying:

- Wipe with dry sponge or soft cloth
- Place the mattress cover in a dry environment (if further drying is needed)
- The mattress cover can be drum dried at temperatures below 50 ° C. When the cover is dry, it must be removed from the dryer immediately, do not leave it in the drum. Interrupting the drying cycle periodically to redistribute the canisters can reduce the drying time (only for bi-elastic canisters).

Disinfection:

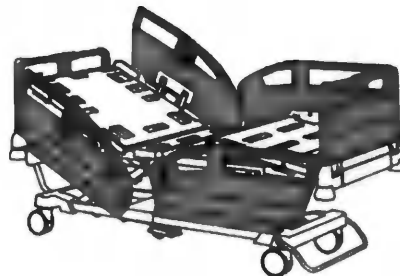
- Wash the surface
- Wipe the surface with a sponge or soft cloth dampened with an approved cleaning agent
- Remove residual water with a soft cloth or sponge
- Dry

General cleaning

We recommend cleaning the bed with detergent and warm water. Do not use an excessive quantity of water, neither abrasive cleaning product.

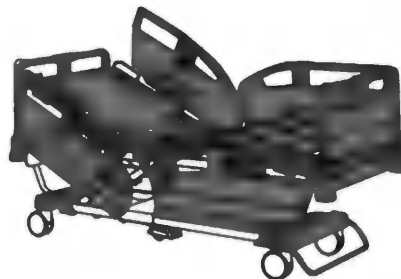
Daily cleaning and disinfection

- Headboard.
- Footboard.
- Rails.

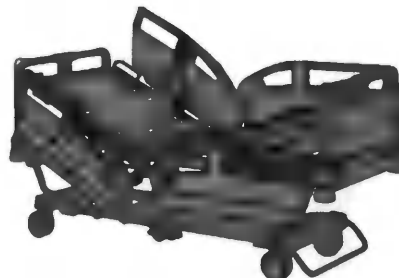
Cleaning and disinfection after a patient leaving

- Headboard.
- Footboard.
- Rails.
- Bed base panels.
- Rails system.
- Protection discs and table.
- Fainting.
- RCP levers.
- Drainages rods.

Exhaustive cleaning and disinfection each two months

- The entire bed.

Steam cleaning

Do not to use any steaming cleaning device to clean NewCare bed range. An excessive humidity may damage the bed mechanisms and components.

Stubborn stains cleaning

To clean stubborn or indelible stains, we recommend using standard domestic cleaning products and a soft-bristled brush.

To clean thick or dry stains, it may be necessary to saturate them first.

Disinfection

Dilute the disinfectant and/or germicides according to the instructions on the manufacturer's label.

Apply the product with a rag avoiding an excessive amount of liquid, or use spray on the surface to be disinfected and spread the product with a rag.

Let it dry for the time recommended by the manufacturer of the product to guarantee maximum effectiveness. We recommend keeping track of the disinfections carried out on each bed, collecting, in a document or record the date of completion, the identification of the bed, the time elapsed since the last disinfection and the procedure and / or products used.

Or if none of the approved cleaning agents provide the required level of disinfection, use a 5.25% sodium hypochlorite solution diluted as follows: 1 part disinfectant to 100 parts water (by volume).

Failure to follow the instructions above may void the hardware warranty.

Bed maintenance

Bed storing

For a long storing period, the bed has to be in its lower height, covered, vibrations free and at a temperature between 0° and 40° C.

Do not use the bed as support (for working or storing purposes) throughout the storage period.

For the correct battery functioning, it is recommended to charge it every two months whilst the storing period.

After a storage period, it is recommended to plug the bed for 10 hours for a total battery charge before putting the bed into service.

Let the bed reach the room temperature before using it after a long storing period.

Preventive maintenance

The bed should be submitted to an efficient maintenance program. This will ensure the operating life and product productivity, and it will help to reduce the deactivation time due to an excessive use.

Every two years, a qualified client service technician should check the electrical equipment to ensure the correct functioning condition.

Preventive maintenance guide

The preventive maintenance guide serves to help the technician to do the correct bed maintenance recommended procedures.

All the points have to be followed making the appropriate and necessary changes during the procedures. This guide is designed to be used along with the preventive maintenance chart. This chart is designed to keep a maintenance history and the consequent repair costs of the bed. Taking specific notes and fulfilling the maintenance can increase the bed life time.

Bed overall status (check every two years)

- Check the general status of the bed verifying there are no visible structural defects such as dents or alarming deformations.
- Verify there are no corrosion signs with extra care on welding joints.
- Observe that all pieces are correctly fitted.
- Check the identification labels are in good conditions and well placed; replace them if necessary.

Make all the repairs, paintwork or component replacement in case it is necessary.

Headboard and footboard (check every two years)

- Check the correct status of the headboard and footboard and its fixations.
- Ensure the locking system of the bed is in good use conditions and locks the footboard correctly.

Repair all the things and replace the components if necessary.

Rails (check every two years)

- Verify that the rails are in good condition and there are not kinked or bent.
- Check that stays properly locked in its highest position.
- Check if the rail releasing system is in good state of use.
- Check that, lowering the rail, the absorber operates correctly.

Make all the needed repairs and replace the components if necessary.

Emergency RCP (check every two years)

- Lift electrically the backrest until the maximum tilted position is reached.
- Verify that the cables that goes from the emergency RCP actioning levers to the motor are in good condition.
- Action independently one lever and check the correct lowering of the backrest.
- Perform steps 1 and 4 activating in step 4 the other lever.

Make all the needed repairs and replace the components if necessary.

Universal positioning of the feet module (check every two years)

- With the electric control rise the legs module until the maximal tilting point is reached.
- Verify the correct condition of the absorbers and the proper fixation with the bed.
- Lift the feet module to test the correct operation.

Make all the needed repairs and replace the components if necessary.

Wheels (check every two years)

- Check the condition of each one of the wheels and that rotates freely.
- Verify the correct condition of the tread and that it has not erosion signs.

Make all the needed repairs and replace the components if necessary.

Multifunction pedal (check every two years)

- Step on the multifunction pedal until it reaches its lowest position (beds brake position).
- Verify that the bed does not move from its place and all the wheels are locked.
- Rise the multifunction pedal until its middle position (free movement position).
- Move the bed and verify that it can make the movements and that no wheel remain locked.
- Rise the multifunction pedal to its highest position (directional fifth wheel position).
- Verify that the fifth wheel is active (it has to touch the floor without forcing the mechanism) and move the bed to ensure it turns on itself without lateral movements.

Make all the needed repairs and replace the components if necessary.

Power supply cable (check every year)

- Verify the cable is not torn or cut.

Make all the needed repairs and replace the components if necessary.

Patient wired control (check every two years)

- Visually check the correct condition of the control plug.
- Verify the wire is in good condition.
- Confirm the correct buttons status.
- Select for two seconds each one of the movements that can be done with the control.
- It must perform the selected function and it has to be a continuous motion.

Make all the needed repairs and replace the components if necessary.

Membrane controls (check every two years)

- Visually check the correct condition of the control.
- Check the membrane is not worn and it is well attached.
- Verify the correct condition of the control plug.
- Select for two seconds each one of the movements that can be done with the control. It must perform the selected function and it has to be a continuous motion.

Make all the needed repairs and replace the components if necessary.

Nurse control (check every two years)

- Visually check the correct condition of the control.
- Check the membrane is not worn and it is well attached.
- Verify the correct condition of the control plug.
- Select for two seconds each one of the movements that can be done with the control. It must perform the selected function and it has to be a continuous motion.
- Check that the control can lock and unlock correctly all the movements.

Make all the needed repairs and replace the components if necessary.

Elevation system (check every six months)

- Check the motor set to ensure if the axis and fixations are correctly assembled.
- Make a whole bed lifting / lowering cycle.
- Verify that is not heard any friction strange noise and there is not an overload signal whilst moving.
- Replace the motor in case of dysfunction.

Make all the needed repairs and replace the components if necessary.

Backrest motor (check every two years)

- Verify the correct motor attachment to the bed.
- Lift the backrest to its highest position and afterward down it to its lowest position (flat).
- Check that no strange friction noise can be

heard and no overload signal sounds whilst moving.

Make all the needed repairs and replace the components if necessary.

Legs motor (check every two years)

- Verify the correct motor attachment to the bed.
- Lift the legs section to its highest position and afterward down it to its lowest position (flat).
- Check that no strange friction noise can be

heard and no overload signal sounds whilst moving.

Make all the needed repairs and replace the components if necessary.

Battery (check every two years)

- If necessary recharge the battery.
- Unplug the bed from mains.
- Try to perform all the beds functions and check that works correctly and there is not disfunctionsigns.
- Check that the battery date is not over three or four years. If it is exceeded, it is recommended to be replaced.

Make all the needed repairs and replace the components if necessary.

Unbrake sound alarm (check every two years)

- Connect the bed to mains.
- Unbrake to hear the alarm.
- Brake again to check the alarm stops.
- Ensure the cable is in correct condition.

Make all the needed repairs and replace the components if necessary.

Rapid problems resolution

PROBLEM	SOLUTION
The bed beeps each time a button is pressed.	It is possible the power is not reaching the bed and it beeps to warn that the battery is discharging.
The bed has locnctions.	Some or all bed functions are locked. To unlock follow the next instructions: From patient control: insert the lock key in the locked function and turn it counter-clockwise until the indicator light (L) turns to green colour. From the nursing control: The unlocking of the bed base sections or bed functions is made by pressing the locking button with a lock icon (K) and following press once the button of the desired section (B, D, H, F, J). The light indicator of the selected section (L) will turn off. From the membrane nursing control on footboard: The unlocking of the bed base sections or bed functions is made by pressing the locking button with a lock icon (K) and following press once the button of the desired section (R, S, U). The light indicator of the selected section (L) will turn off.
Leds blinks and the bed does not work.	Wait till the leds stop blinking and check the locking situation as the preview section.
The bed is locked and nothing works.	Unplug from mains and battery. Wait 30 sec. and plug to mains again. Once the controls are recognised plug the battery.
The bed stay locked and the above procedure did not work.	In this case you need to reset the system, follow these steps: From the patient control: Press and hold A and C keys for a few seconds until the control beep disappears. The system resets and operates correctly again. From the nursing control: Press and hold simultaneously C and D keys for a few seconds till the system resets. From the membrane nursing control on footboard: Press and hold simultaneously 1 and 2 keys until the control beep disappears.
After reset the system the bed operates but when pressing any button it locks again.	If after reset the bed by pressing the controls it locks again, it is possible that there is a problem with the connexion cables and produces a short circuit in the control box and stops the bed. Replace the cable and call the technical service.
The battery does not work.	If after the battery charge, the charge indicator located in the nursing control does not shows that it is charged, check the battery condition. If it is below 17 volts, replace it for a new one.

To avoid the battery discharge, the bed will be shipped with the battery unplugged. That's why, once the bed has been plugged to mains and after checking the controls (nursing control leds stop working), the battery has to be connected to the control box.

The batteries have a limited life that depends on the use done of them. The value from which it is impossible to recharge is 17 volts. Below this value, the battery is unrecoverable and it has to be replaced. Below 20 volts, it is already a low value and it is recommended

The batteries have to stay in the bed when it is plugged to mains. In case the bed is storages, it is recommended to unplug the battery and if the storing time is over 2-3 months, is better to remove the battery and leave it on a charger or use it in another bed that is in plugged, to avoid discharges.

not surpass it.

The beds clutch is only to lower de bed and not to lift it. The inappropriate use of the system voids the backrest motorand clutch warranty.

Once the membranes are stuck, do not reuse them because they can be damaged whilst taking them off. When a membrane is removed, it is important to clean the supporting surface to avoid damages on the pulsation.

The connexion box has an overuse protector (aprox. for 2 minutes) that is added to the reversible and not reversible thermics (80 y 95°) that has the control box as safety elements.

If the connectors are hard to plug, do not remove the o-rings (it may lose the IP66 protection). It should apply a little bit of grease or vaseline to facilitate its connection.

When disconnecting from mains always wait a couple of minutes to connect it again and press the controls buttons todischarge the control PCB capacitors.

Product disposal

NewCare and electronic accessories shouldn't be litter with other urban wastes. To protect natural resour- es and promote recycling, separate batteries and recycle them by local collection procedure. Other elements in NewCare can be separated: all metallic elements, plastic elements, wooden elements, phenolic elements, packaging and RAEES, to avoid environmental damages classify these elements and recycle them properly or entrust them in a residues delivery point in your city. Thus material sources recycling is promoted.

Private users can contact the establishment where they purchased NewCare or competent authority to be informed about how and where they can deliver the bed to be recycled safely. Please comply with the EU Waste Electrical and Electronic Equipment Directive.

All electronic equipment complies with the Restriction of the Use of Certain Hazardous Substances in Electrical and Electronic Equipment Directive.

The medical device and / or its parts must be disposed of in accordance with the existing applicable local regulations.

The bed must be cleaned and disinfected before disposal.

~~The bed and / or its parts and accessories must be disposed of in accordance with the existing~~

applicable local regulations. For the Russian Federation, it is recommended to dispose of in accordance with the rules approved by SanPiN 2.1.3684-21 Sanitary and epidemiological requirements for the maintenance of the territories of urban and rural settlements, for water bodies, drinking water and drinking water supply, atmospheric air, soils, living quarters, the operation of industrial, public premises, the organization and implementation of sanitary and anti-epidemic (preventive) measures for class "A »(Non-hazardous) waste, excluding electronic parts.

Electronic parts must be disposed of in accordance with local special regulations

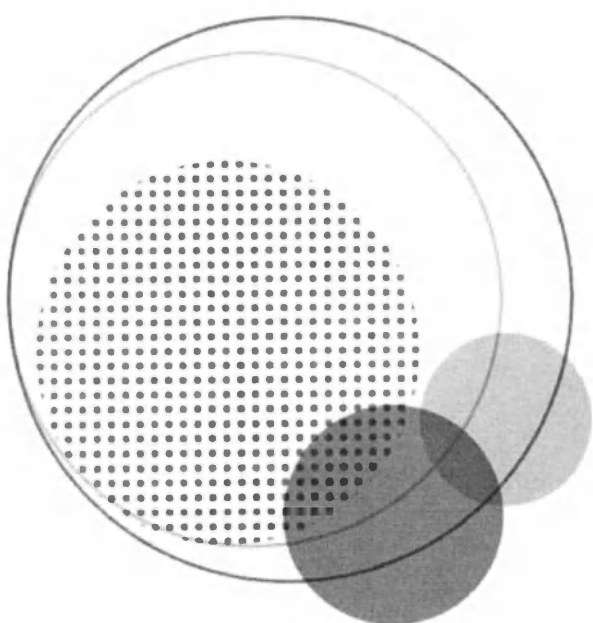
This symbol indicates a separate collection of waste for electrical and electronic equipment

Warranty

I.H. Pardo will only be responsible for the safety and reliability of those products that have been periodically re-viewed and used in accordance with safety regulations.

This product is guaranteed for 24 months from the date of purchase. This warranty covers all damage and errors related to materials and workmanship. Damage and errors caused by incorrect use or external causes are not covered by the guarantee.

Technical specifications

According the model bed, the dimensions can be vary.

Physical bed description

Measurement and weight specification	Dimension
Total length with headboard and footboard (extended until 2300 mm useful surface)	
• Standard Model:	2230 mm. $\pm$ 10 mm.
• S Model:	2210 mm. $\pm$ 10 mm.
• XS Model:	2080 mm. $\pm$ 10 mm.
• XSS Model:	2060 mm. $\pm$ 10 mm.
Total width with side rails	
• Standard Model:	980 mm. $\pm$ 10 mm.
• S Model:	960 mm. $\pm$ 10 mm.
• XS Model:	980 mm. $\pm$ 10 mm.
• XSS Model:	960 mm. $\pm$ 10 mm.
Minimum / maximum height of bed base (without mattress) (low - raise position)	
• Integral wheel $\varnothing$ 150 mm:	410 mm. - 815 mm. $\pm$ 10 mm.
• Double wheel $\varnothing$ 150 mm:	400 mm. - 805 mm. $\pm$ 10 mm.
• Metallic wheel $\varnothing$ 150 mm:	415 mm. - 820 mm. $\pm$ 10 mm.
Crane passing clearance height	
• Integral wheel $\varnothing$ 150 mm:	170 mm. $\pm$ 10 mm.
• Double wheel $\varnothing$ 150 mm:	160 mm. $\pm$ 10 mm.
• Metallic wheel $\varnothing$ 150 mm:	175 mm. $\pm$ 10 mm.
Length of extended bed base	300 mm.
Bed base double-regression	150 mm. $\pm$ 10 mm.
Length of bed base sections	
Standard Model:	790 / 170 / 320 / 620 mm.
S Model:	790 / 170 / 320 / 620 mm.
XS Model:	790 / 170 / 320 / 470 mm.
XSS Model:	790 / 170 / 320 / 470 mm.
Angle of backrest tilting	65° $\pm$ 5°
Angle of legrest tilting	45° $\pm$ 2°
Angle of footrest tilting	18° / 7° $\pm$ 2°
Angle of Trend / Reverse trend tilting	13° / 15°
Unladen bed weight	From 125 kg. to 150 kg.
Safe working load (SWL)	250 kg.
Maximum patient weight (MPW)	185 kg.

Rails physical description

Specifications	Dimension
Length of backrest rail	1010 mm.
Length of legrest rail	1040 mm.
Length of legrest rail (safe exit)	719 mm.
Rails height respecting the bed base	430 mm.

Recommended mattress description

Specifications	Dimension
Width x Length of mattress (recommended)	
• Standard Model:	850 x 2000 mm.
• S Model:	800 x 2000 mm.
• XS Model:	850 x 1900 mm.
• XSS Model:	800 x 1900 mm.
Minimum height	130 mm.
Maximum height	210 mm.

Use environment (according to direction EN 60601-2-52:2010)

Use environment 1	Intensive care hospital.
Use environment 2	Acute care hospital.
Use environment 3	Long stay hospitalization.
Use environment 5	Outpatient care.

Use environment physical description

Condition	Range
Temperature	Between +5°C and +40°C
Humidity	
Pressure	From 500 hPa to 1060 hPa (from 0 to 2000 m altitude)

Storage environment physical description

Condition	Range
Temperature	Between +5°C and +40°C
Humidity	20% - 80% at 30°C
Pressure	From 500 hPa to 1060 hPa (from 0 to 2000 m. altitude)

Electric system description

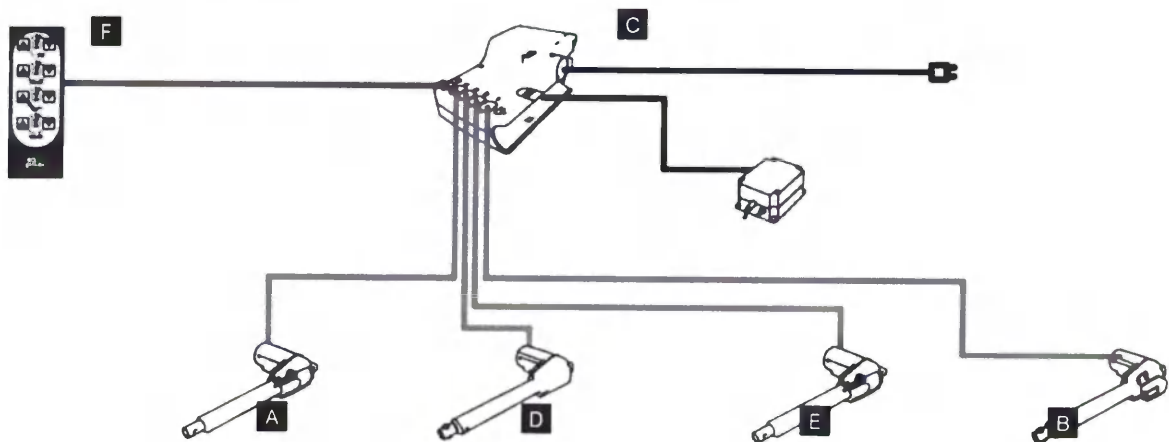
Description	Specification
Electric shocks protection type (CEI 60601-1)	Class II
Electric shocks protection grade (CEI 60601-1)	Type B
Water entry protection (CEI 60529)	IPX6

Electrical diagram

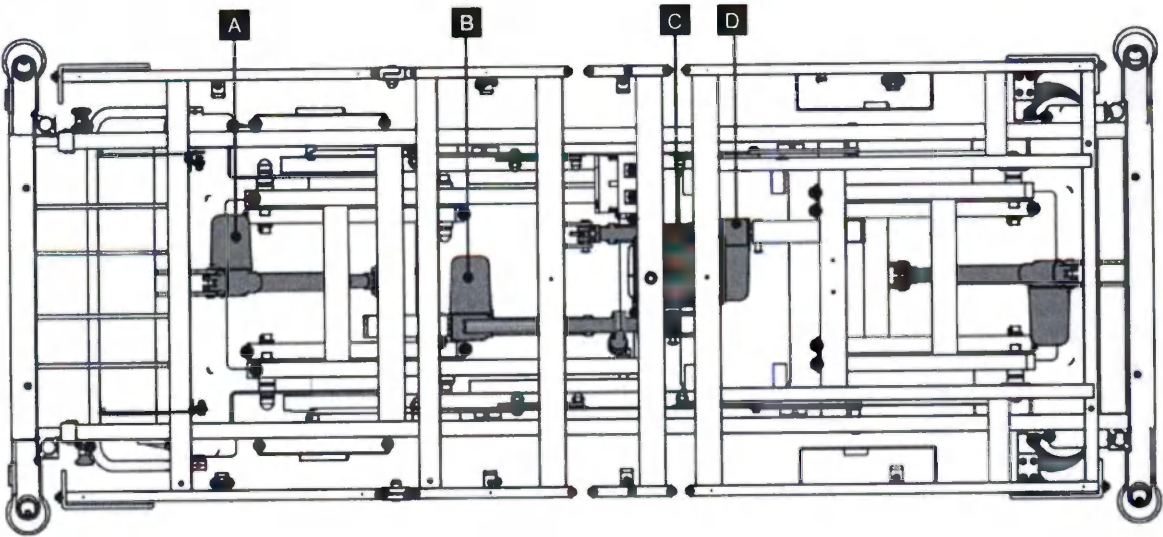
NewCare bed has the following basic configuration:

- Control box.
- Battery.
- Four linear motors one for the backrest, one for the legs and two for the lifting).
- Patient control with extendible wire (depending on options).

The wiring diagram* is the following:

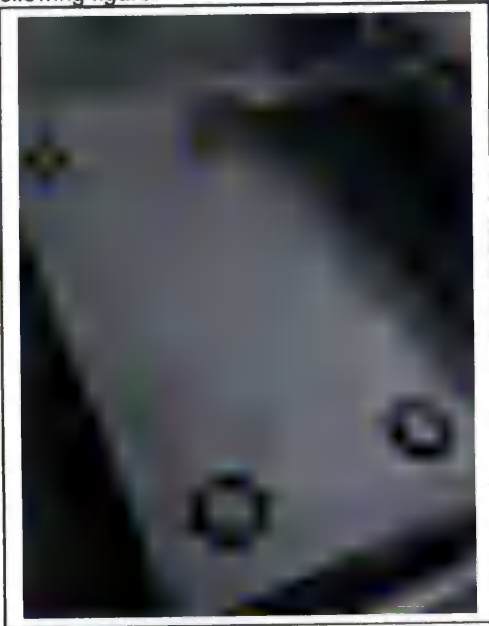
*This diagram is just an example of the possible configuration.
The complete diagrams of each bed configuration are in the spare parts manual.

Item	Description	Item	Description
A	Bed lifting motor (legrest section)	D	Legrest lifting motor
B	Backrest lifting motor	E	Bed lifting motor (backrest section)
C	Control box	F	Patient control

Packaging of the bed

- If specified in the contract or order, the medical beds will be packed according to the model called "double lid boxes".
- The pallet is placed on the floor, one of the covers is placed on top, after closing the sashes are stapled or glued.
- Depending on the size of the bed, the mooring (fixing) wheel cylinders are placed on the first cover, as shown in the following figure:

- The bed is positioned to align with the cylinders as follows:

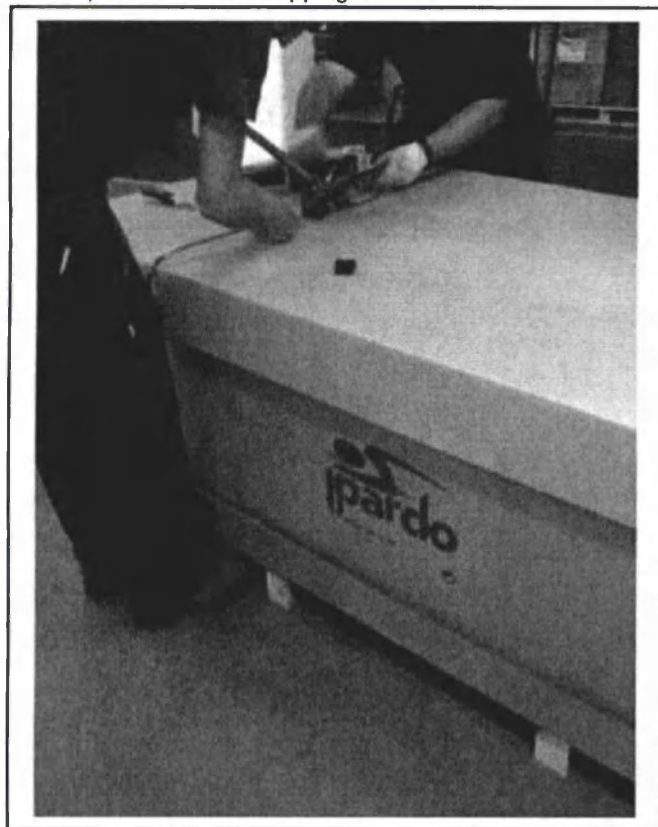
- The cylinders are screwed to the base.
- The side of the box is positioned as follows:

- The cover is fitted the same size as the base

- Then proceed to the strapping and identification of the kit.

Installation of the bed

- By definition, installation is the last stage of the sale. It is carried out at the client's premises and mainly with the use of material resources and human resources, which are managed by the company Industrias Hidraulikas Pardo C.L.

- For this reason, the technical service plans the specified installation as follows: clarify with the customer all data that may be of interest for the installation (for example, the specific day and time of delivery, truck access to the Center, documentation requested by the customer, etc.) . All this information is agreed between the customer and the installer.

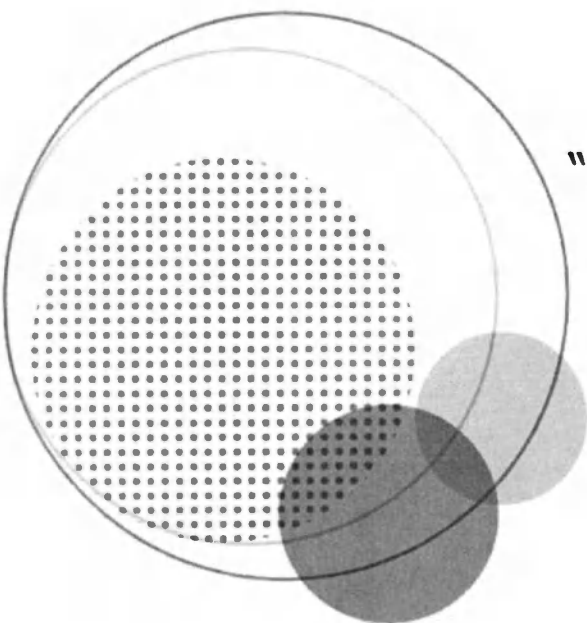
All additional information, where appropriate, is attached to the rest of the documentation that will be delivered to those assigned to this external assembly. The specified documentation will also contain at least "waybills". 3 copies of invoices are delivered: one - to the customer as acceptance of delivery, another signed by the driver as confirmation of delivery, and the third is returned by the installer as confirmation of installation.

If the exterior installation is performed by the service provider, then it must be properly assessed and approved as such. In the case of using the services of temporary employment agencies, the staff of the Pardo company will be present (commercial, quality, SAT, etc.).

- The person doing the programming owns his own funds and responsible employees, and freely contracts with approved suppliers.

(Carriers, temporary work companies, etc., aids they deem suitable for the installation).

- In addition, installation means that the product is definitely suitable for the customer's use.

"EC" Declaration of conformity

"EC" DECLARATION OF CONFORMITY

INDUSTRIAS HIDRAÚLICAS PARDO, S.L.

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50198 La Muela, Zaragoza (Spain) CIF:

B50057322

Authorized company by the Spanish Medical and Sanitary Products Agency of the Spanish Health Ministry
with num. 261 P.S.

STATES THAT THE PRODUCTS DESIGN AND MANUFACTURING ARE:

Type:	Hospital bed		
Family:	NewCare	Basic UDI-DI	8431212CAMA-HOSPITAL-NCAG
Model:	NewCare, NewCare S, NewCare XS, NewCare XSS		
Serial Num.:	See Labels	Manufacturing date	See Labels
Classification:	These products are class I, according to Chapter V of Regulation (UE) 2017/745		

IN ACCORDANCE WITH THE FOLLOWING DISPOSITIONS:

- Regulation 2017/745 of Parliament and the Council, of 05 April 2017, concerning medical devices.
- Directive 2006/42/CE of Parliament and the Council on 17 May 2006 concerning machines and amending the Directive 95/16/CE.
- Directive 2006/95/CE of Parliament and the Council, of 12 December 2006 concerning on the approximation of the laws of the Member States relating to the electrical material to be use don certain voltage limits.
- Directive 2004/108/CE of Parliament and the Council, 15 December 2004 concerning on the approximation of the laws of the Member States relating to the electromagnetic compatibility.

REFERENCE STANDARDS:

- Standard UNE EN ISO 12100, "Machinery safety. General principles for design, risks analysis and risks reduction".
- Standard UNE 20324, "Degrees of protection provided by enclosures (IP Code)". Equate to IEC 60529"
- Standard UNE EN 60204-1, "Machinery safety. Electric equipment of the machines. Part 1: General requirements"
- Standard UNE EN 60601-1, "Medical electrical equipment. Part 1: Basic safety general requirements and essential functioning"
- Standard UNE EN 60601-1-2, "Medical electrical equipment. Part 2: Safety general requirements. Collateral standard: Electromagnetic compatibility. Requirements and tests"
- Standard UNE EN 60601-2-52, Medical electrical equipment. Part 2-52: Basic safety general requirements and essential functioning of hospital beds"

Adolfo Gaudes
TECHNICAL QUALITY MANAGER
Zaragoza, 19th March 2020

Client:

Model:

Serial

KEY	
L = Lubrication	R= Repair Replacement
M = Cleaning	Ok = No problem
A = Adjustment	NA = Not applicable

num.:

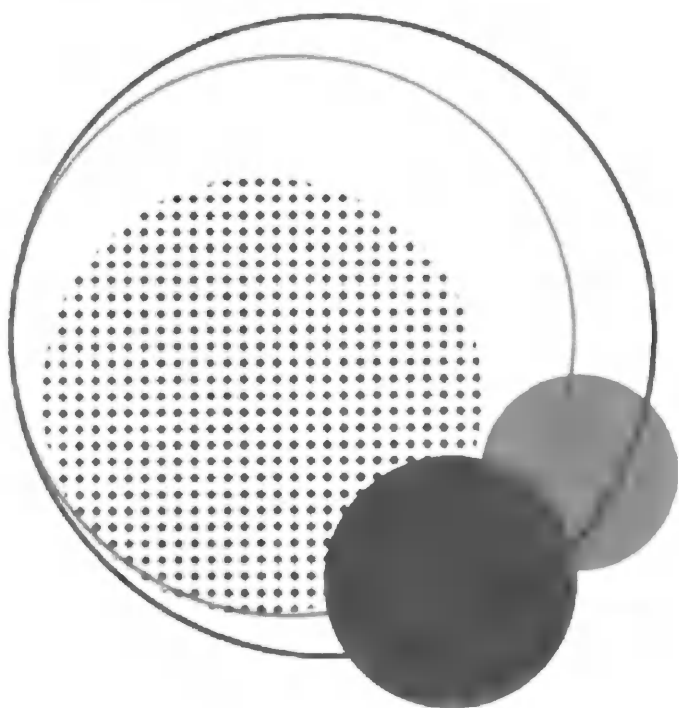
FUNCTIONS								
	Date	Date	Date	Date	Date	Date	Date	Date
General Status								
Headboard and footboard								
Brake system								
Electric system								
Rails								
Wheels								
Bed base								
Sections motors								
Lifting motors								
Extensible, Spring gas, CPR								
Others								

Repair cost								
-------------	--	--	--	--	--	--	--	--

Labour cost								
-------------	--	--	--	--	--	--	--	--

Total cost								
------------	--	--	--	--	--	--	--	--

Reviewed by								
-------------	--	--	--	--	--	--	--	--

Matresses

All bed models can be equipped with the following mattress models:

Model	Q-ty of layers	Foam type	Уровень защиты от пролежней	Max patient weight	Top profile	Fire resistance	Symmetry
Density 900	3	HR + Visco	I – very high	250 kg	Yes	No	Yes
Density 650 Fresh	3	HR + Visco	II – high	250 kg	No	Yes	No
Density 600	3	HR + Visco	II – high	250 kg	Yes	No	Yes
Density 475 Fresh	2	HR + Visco	II – high	230 kg	Yes	No	Yes
Density 450	2	HR + Visco	III-middle	230 kg	Yes	Yes	Yes
Density 400	2	HR + Visco	III-middle	230 kg	No	No	Yes
Density 250 Fresh	1	Biocell	IV – low	190 kg	No	No	Yes
Density 150	1	HR	IV – low	190 kg	Yes	No	Yes
Density 100	1	HR	IV – low	190 kg	Yes	No	Yes

HR - Polyurethane foam (latex free)

Visco - Memory Foam (Latex Free)

Biocell - Biocell foam BIOCEL3230 (latex free) Eco-friendly foam with high elasticity, made from natural plant products (in particular soy), renewable and environmentally friendly.

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Colchón hospitalario / Hospital mattress

DENSITY 100

FIRA
Firearm Safety Research Association

CE

Colchón hospitalario / Hospital mattress

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Mattress specially designed for ambulatory or hospital patients with low ulceration risk, as well as for the treatment of early stage ulcers and sores.

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High density monolithic foam core (25 kg/m^3), with open cell structure that improve the air exhaust through the foam, becoming the core breathable and hygienic.

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14 cm. core height, which avoids the contact between patient and bed, therefore, preventing the pressure contacts.

The core structure and composition makes this mattress suitable for all sanitary areas.

Fireproof core according to the "Furniture and Furnishings (Fire)", BS 5852:2 CRIB2, UNE EN 597-1 and EN 597-2 regulations.

Ergonomic mattress, keeps the body in his natural curved position, relieved, preventing the back and cervical pain.

Reversible and symmetric core allows using indistinctly the head or feet areas.

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This mattress is not recommended for patients considered to be bariatric.

It is indicated for patients with risk of LOW ulceration.

This mattress is designed for both beds with phenolic beds and beds with beds of slats.

Use each bed mattress with the appropriate bed measurement. Do not use on smaller beds.

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Ulceration risk*

Bajo

*According to Norton equals degree of ulceration II

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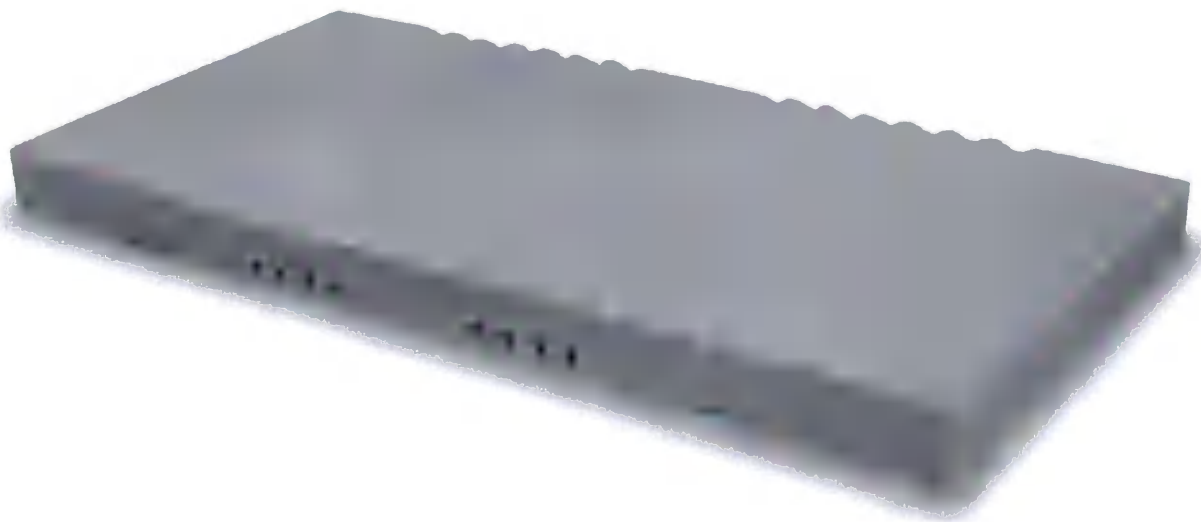
Specially designed mattress to be used in hospital and health centres, for patients with low ulceration risk or for mitigate up to grade I pressure ulcers or sores.

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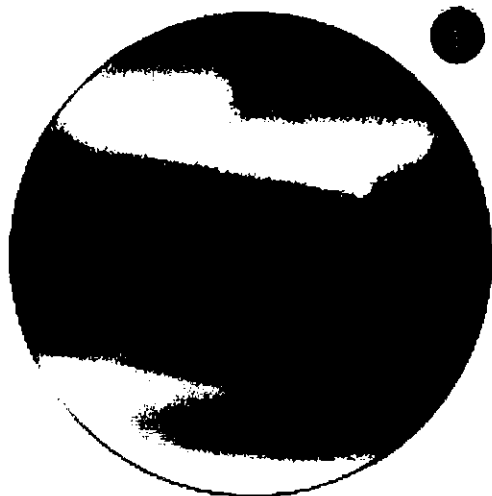
High resilience polyurethane foam (30 kg/m³), which open cells structure, AIRFLOW system and mattress perforations, provide a high lightness; 14 cm. thick.

High performance AIRFLOW system, it allows a great adaptability and breathability.

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Corrugated upper layer, surface divided in 5 sections. The corrugated areas are placed in the shoulders and hips sections, to improve the progressive adaptability of the patient.

Framed inferior layer, in 5 sections, for better adaptability and movement of the sections

Fireproof core according to the "Furniture and Furnishings (Fire)", BS 5852:2 CRIB2, UNE EN 597-1 and EN 597-2 regulations.

All foams are Oeko-Tex class I and REACH certified.

Symmetric mattress allows using indistinctly the head or feet areas.

OPTIONAL: Cut the mattress corners for better adaptability to beds.

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Specially designed mattress to be used in hospital and health centres, for patients with low ulceration risk and for the treatment of early stage ulcers or sores.

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High resiliency ecologic Biocell foam (32 kg/m³), made of vegetal origin natural products (specifically soya), renewable and non-polluting; with a height of 15 cm. thick.

Core is composed by open micro cells, obtained from open seeds, ensuring a maximal breathability.

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Framed lower mattress layer for an optimal patient pressure whilst his bed placement. High adaptability of the mattress to the bed, promoting its articulation.

Totally ergonomic, keeps the body in his natural curved position, relieved, preventing the back and cervical pain.

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Symmetric core allows using indistinctly the head or feet areas

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Ulceration risk*
Bajo
*According to Norton equals degree of ulceration II

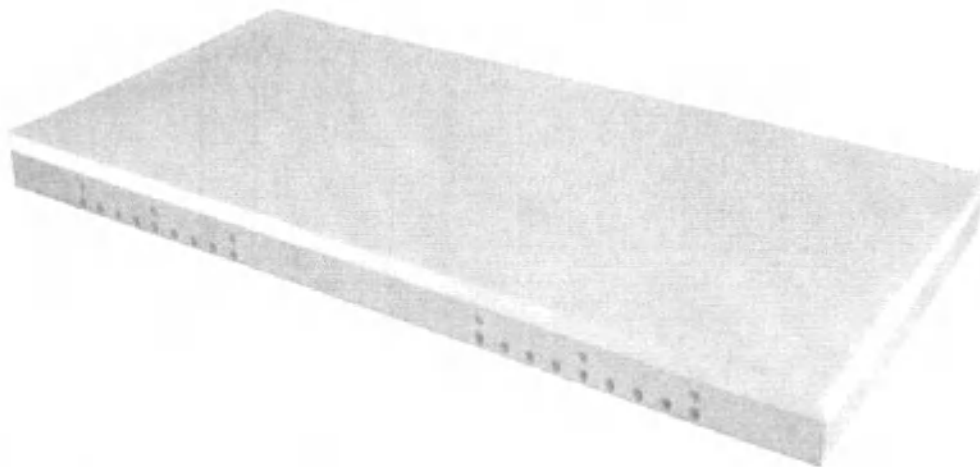
Espuma Biocell	Densidad Biocell	Libre de látex natural	Origen vegetal	Material ergonómico	Material transpirable	Material antialérgico
Biocell Foam	Biocell density	Natural latex free	Vegetal origin	Ergonomic material	Breathable material	Anti-allergy material

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Specially designed mattress to be used in hospital and health centres, for patients with medium-low ulceration risk or for patients with medium grade ulcers or sores.

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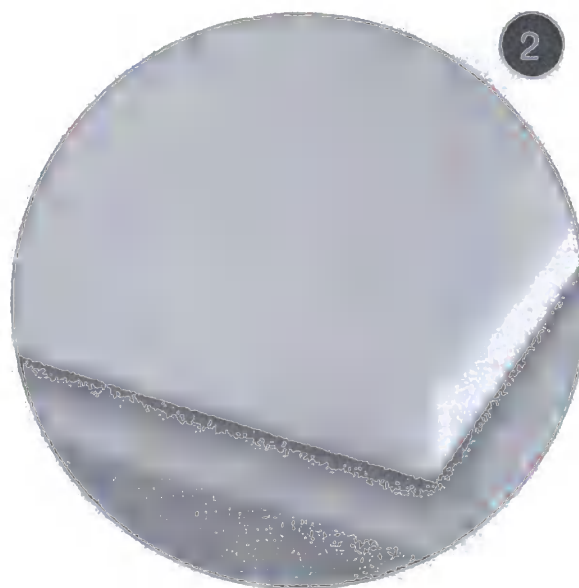
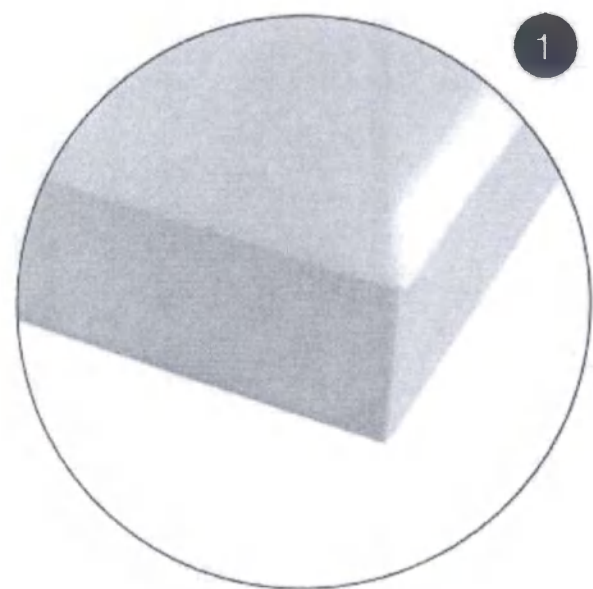
Core with two layers, base with HR polyurethane foam base freefoam (28 kg/m³). 9 cm. height. Due to its porosity, allows a high breathability and hygiene.

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Upper layer made of high density thermo sensible viscoelastic foam 50 kg/m³; 5 cm. height. The Viscoelastic foam increases the patient adaptation capacity reaching low capillary occlusion levels, avoiding the PU's appearance.

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Framed lower mattress layer for an optimal patient pressure against the bed placement. High adaptability of the mattress to the bed, promoting its articulation.

Suitable for patients with medium-low ulceration risk and not considered bariatric.

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It is indicated for patients with risk of ulceration MEDIUM - LOW.

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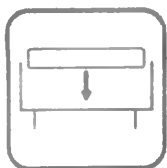
Specially designed mattress to be used in hospital and health centres, for patients with medium-high ulceration risk, as well as for mitigate pressure ulcers or sores in medium stage.

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High density polyurethane, Biogreen, foam based core (31 kg/m^3). Made of vegetal origin natural products, renewable and non-polluting. 10 cm. height.

Biogreen is composed by open micro cells, which provide maximum breathability.

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Upper layer made of high density (50 kg/m^3) MF viscoelastic. Bio, which combines the Biogreen foam properties and the MF foam. Mind Foam (slow memory and high density viscoelastic foam, with a wrap around effect). 5 cm height.

Enhanced core with Aloe Vera extracts that provide freshness.

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Upper layer conformed with different profiles, differentiating 5 sections, according to the patient's pressure areas, to improve the blood circulation and pressure of the patient.

Mattress base framed, for optimal bed core adaptability, when the bed sections are articulated.

Microperforations made in the viscoelastic layer, to provide greater mattress breathability.

Comfortable, ergonomic and non-deformable

All foams are Oeko-Tex class I and REACH certified

Symmetric core allows using indistinctly the head or feet area

OPTIONAL: Cut the mattress corners for better adaptability to beds

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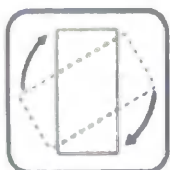
It is indicated for patients with risk of MEDIUM - HIGH ulceration.

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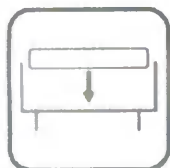
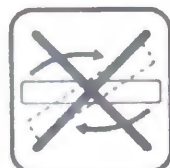
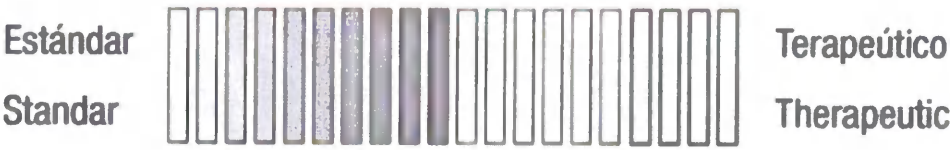
Ulceration risk*
Medium - High

*According to Norton equals degree of ulceration III

							
Poliuretano HR + Viscoelástica	Densidad Espuma HR	Densidad Viscoelástica	Libra de látex natural	Material termo adap- table	Material transpirable	Material antialérgico	Material ignífugo
Polyurethane HR + Viscoe- lastic	HR foam density	Viscoelastic density	Natural latex tree	Thermo- adaptable material	Breathable material	Anti-allergy material	Fireproof material

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Sandwich type mattress for high ulceration risk patients, laying in bed for at least 15 hours per day, as well as to mitigate the pressure ulcers or sores up to grade III.

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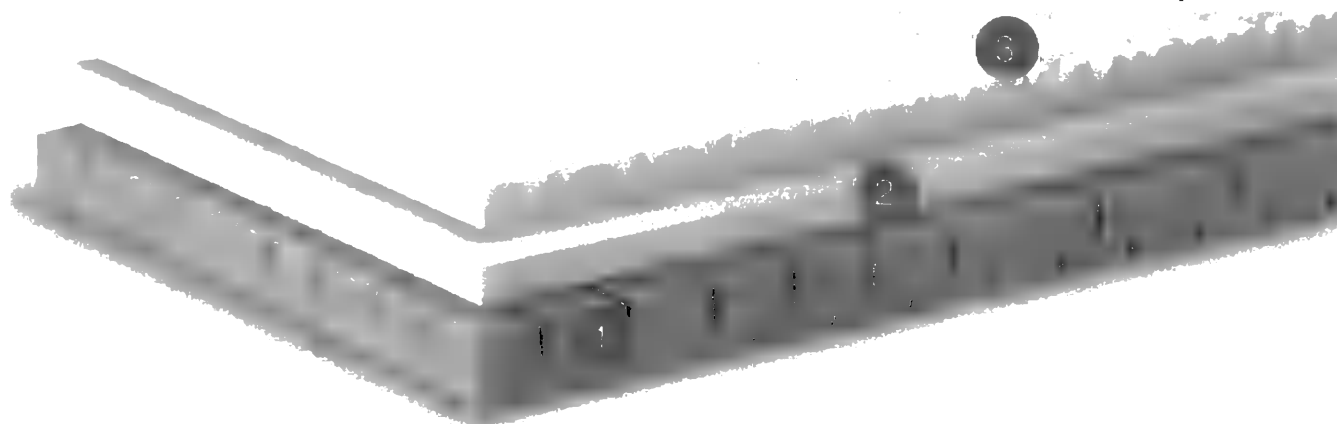
HR Polyurethane foam base (38 kg/m³), with 8 cm. of height. Its height porosity allows a great air circulation through the inner.

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Intermediate layer made of HR polyurethane foam with more density (42 kg/m³) with 4 cm. of height that provides more comfort and stability.

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Thermosensitive viscoelastic foam with 4 cm. height for the upper layer. In contact with the body, the viscoelastic foam turns softer, keeping the same supporting level and ensuring an excellent comfort levels.

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Upper layer conformed with different patterns, divided in 5 sections, depending in patients pressure areas, to improve the blood circulation and pressure of the patient.

Mattress base framed, for an optimal adaptability of the mattress core.

Transversal frame in the upper side of the base, for better patient comfort and core elasticity.

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Sandwich type mattress for high ulceration risk patients, laying in bed for at least 15 hours per day, as well as to palliate pressure ulcers or sores at intermediate stage.

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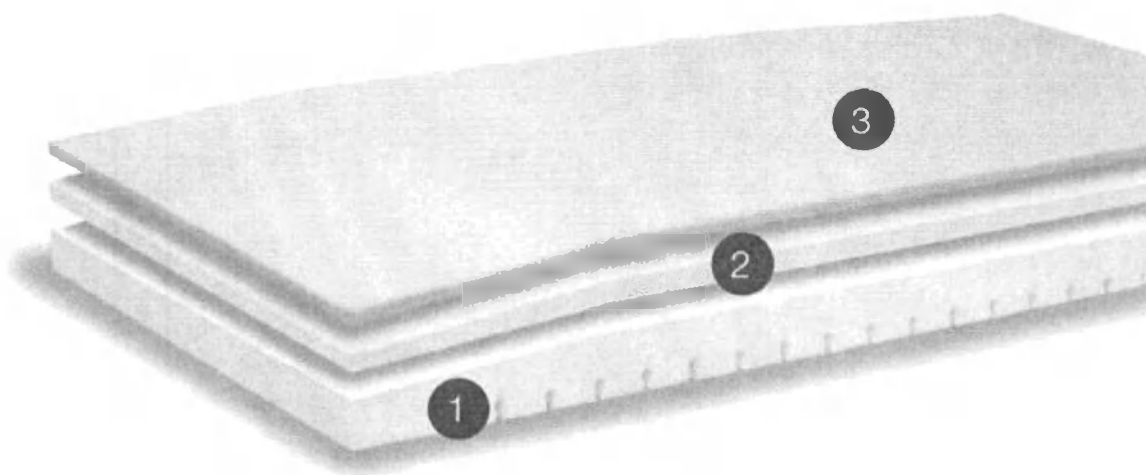
HR polyurethane foam, Freefoam (30 kg/m³), with open cells, that allows high air circulation through the foam, with a great breathability; 10 cm. of height.

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Intermediate layer made of thermofoam viscolastic (50 kg/m³), with 3 cm. height, which improve the adaptation capacity to the patient and his comfort.

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Upper layer in Rafrechi® (50 kg/m³) viscolastic foam, as well as the intermediate layer, but in this case the foam is refreshed for a greater breathability and to avoid patient maceration; 2 cm. height.

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Mattress base frame, for an optimal core adaptability, when the bed sections are articulated.

Non-inflammable core according to the "Furniture and Furnishings (Fire)", BS 5852:2 CRIB2, UNE EN 597-1 y EN 597-2. regulations

Heel-Relief area slightly decreased down area, to prevent the pressure increase this area and avoiding the pressure ulcers generation.

All foams are Oeko-Tex class I y REACH certified.

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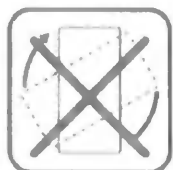
This mattress is recommended for patients not considered to be bariatric.

It is indicated for patients with risk of ulceration HIGH.

This mattress is **designed** for both beds with phenolic beds and beds with beds of slats.

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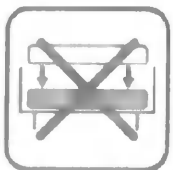
The head and foot area can't be used indistinctly.

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Always place on top of bed.

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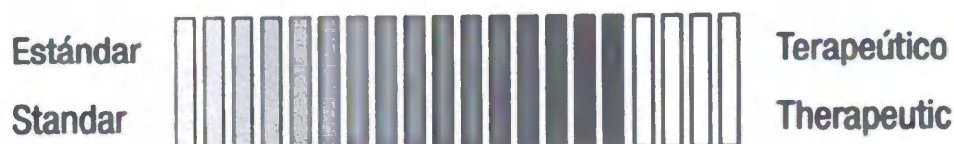
Do not place on top of another mattress.

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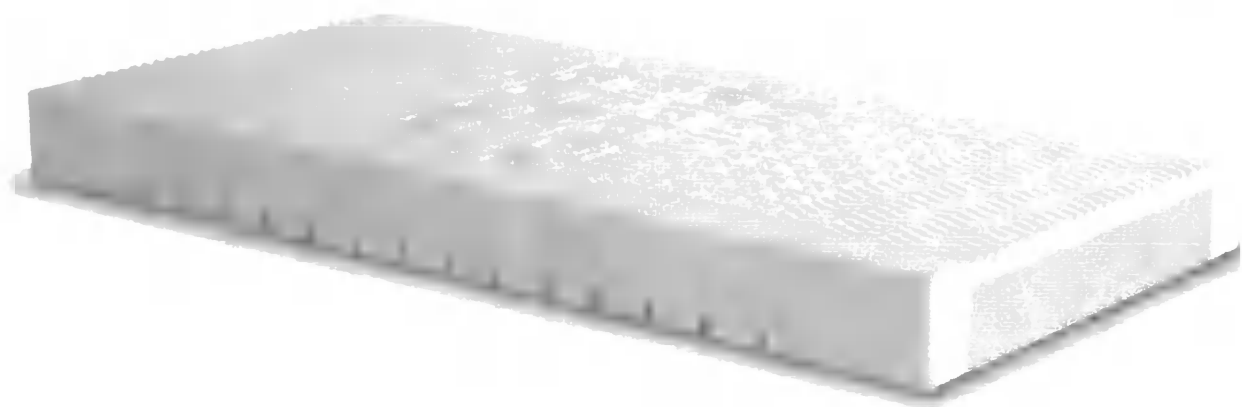
Encapsulated mattress developed for therapeutic use in hospital environment, for patients with chronic pain and very high ulceration risk, as well as for mitigate pressure ulcers or sores up to stage III.

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Polyurethane foam HR (30 kg/m³) lower layer. Base with 8 cm. height. High porosity, which allows a high breathability and hygiene, avoiding the humidity and maceration.

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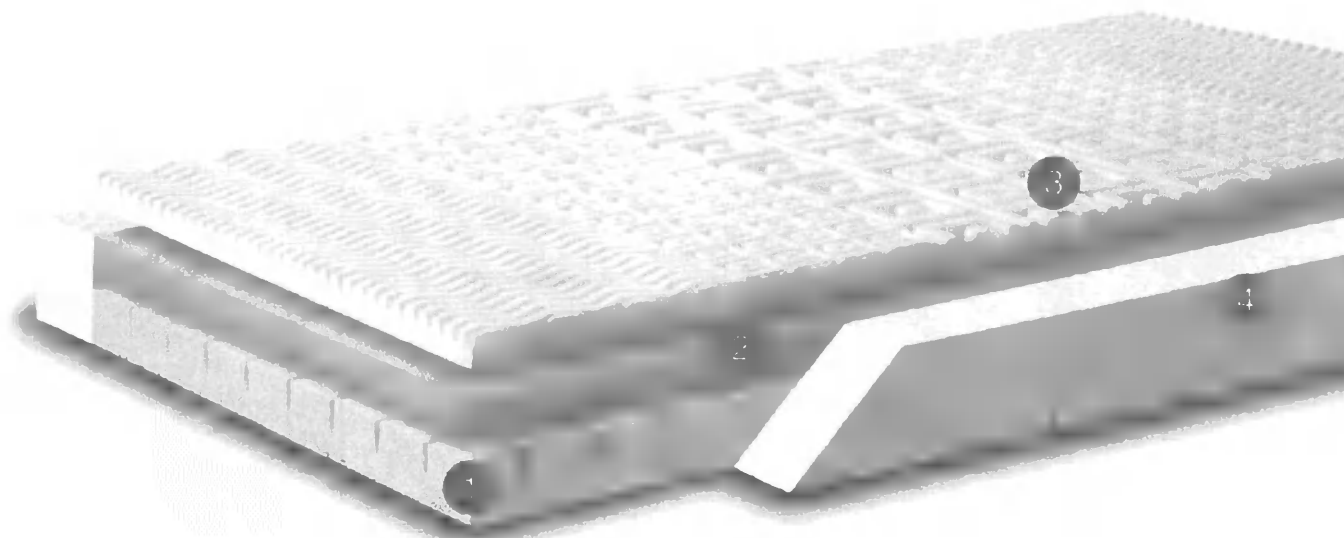
Intermediate layer on HR Polyurethane foam of 4 cm. height. With more density than the base (35 kg/m³), that provides better comfort and stability.

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Upper layer made of thermo sensitive viscolastic foam (50kg./m³), 4 cm. height. In contact with the body, the viscolastic foam turns softer, keeping the same supporting and endowing an excellent comfort levels.

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HR polyurethane foam (40 kg/m³) frame. With 16 cm. height and 10 cm. width.

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The foam frame reinforce the mattress, ensuring the patient stability whilst turning on himself, guaranteeing the durability against the perimeter efforts and improving the comfort when sitting at the bed side.

Upper layer conformed with different patterns, divided in 5 sections, depending the patients pressure areas, to improve the blood circulation and pressure of the patient.

Transversal frame in the upper side of the base, for better comfort, increased the elasticity and breathability of the core.

Mattress base frame, for an optimal adaptability of the mattress core with the bed

All the foams are Oeko-Tex class I and REACH certified.

Symmetric core allows using indistinctly the head or feet areas.

OPTIONAL: Cut the mattress corners, for better adaptability to beds.

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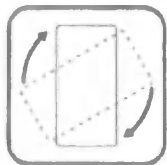
It is indicated for patients with very high risk of ulceration.

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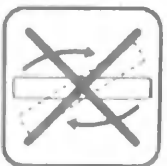
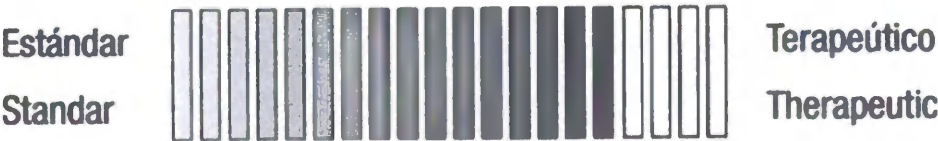
Ulceration risk*
Very High
*According to Norton equals degree of ulceration III

									
Poliuretano HR + Visco Polyurethane HR + Visco	Densidad Inferior Lower density	Densidad Intermedia Medium density	Densidad Carenado Frame density	Densidad Viscoelástica Viscoelastic density	Libre de látex natural Natural latex free	Material termo adaptable Thermo- adaptable material	Material conformado Conformed material	Material transpirable Breathable material	Material antialérgico Anti-allergy material

Model versions and composition

Medical Device: Medical functional electric 4-section NewCare bed in versions with accessories

Versions

1.1. Medical functional electric 4-section NewCare Standard bed consisting of:

1. Bed frame - 1 pc.
2. Bed platform, composition:
 - removable sections of the stock made of HPL-plastic - 3 pcs.
 - fixed section of the stock made of steel - 1 pc.
3. Plastic side rails of the back section (if necessary), execution options:
 - with built-in control panel (if necessary) - 2 pcs.
 - with built-in backlit control panel (if necessary) - 2 pcs.
 - without control panel (if necessary) - 2 pcs.
4. Plastic side rails of the leg section (if necessary) - 2 pcs.
5. Metal side rails, execution options:
 - stainless steel (if necessary) - 2 pcs.
 - painted steel (if necessary) - 2 pcs.
6. Head end, execution options:
 - removable with hand grips (if necessary) - 1pc.
 - removable straight (if necessary) - 1pc.
 - non-removable (if necessary) - 1pc.
7. Foot end, execution options:
 - removable with hand grips with built-in control panel (if necessary) - 1 pc.
 - removable with hand grips without a control panel (if necessary) - 1 pc.
 - removable straight line with built-in control panel (if necessary) - 1pc.
 - removable straight line without control panel (if necessary) - 1pc.
8. Wheels, execution options:
 - metal 150 mm (if necessary) - 4 pcs.
 - integral 150 mm (if necessary) - 4 pcs.
 - double 150 mm (if necessary) - 4 pcs.
9. Central brake, execution options:
 - multifunctional pedal (if necessary) - 1 pc.
 - side pedals of the central brake on the wheels at the foot end (if necessary) - 2 pcs.
 - side pedals of the central brake for all wheels (if necessary) - 4 pcs.
10. 5th guide wheel 100 mm (if necessary) - 1 pc.

11. Foot control pedal (if necessary) - 1 pc

12. Anti-decubitus mattress, variants:

- Density 100 (if necessary) - 1 pc.
- Density 150 (if necessary) - 1 pc.
- Density 250 fresh (if necessary) - 1 pc.
- Density 400 (if necessary) - 1 pc.
- Density 450 (if necessary) - 1 pc.
- Density 475 fresh (if necessary) - 1 pc.
- Density 600 (if necessary) - 1 pc.
- Density 650 fresh (if necessary) - 1 pc.
- Density 900 (if necessary) - 1 pc.

13. Hygienic mattress cover, variants:

- Cugar (if necessary) - 1 pc.
- Select (if necessary) - 1 pc.

14. Additional mattress for extending the bed with a hygienic cover (if necessary) - 1 pc.

15. Lever for emergency lowering of the back section - 2 pcs.

16. Protective bumpers, execution options:

- diameter 84 mm (if necessary) - 4 pcs.
- diameter 104 mm (if necessary) - 4 pcs.

17. X-ray cassette holder with a guide (if necessary) - 1 pc.

18. Electric drive - 4 pcs.

19. Pneumatic leg section adjustment - 1 pc.

20. Mattress holders - not more than 4 pcs.

21. Mounting holes - not more than 4 pcs.

22. Control panel for medical personnel (if necessary) - 1 pc.

23. Patient control panel, execution options:

- backlit, without blocking of functions (if necessary) - 1 pc.
- without blocking of functions (if necessary) - 1 pc.
- with blocking of functions (if necessary) - 1 pc.

24. Power cable - 1 pc.

25. Battery (if necessary) - 1 pc.

26. Built-in segment for lengthening the bed 300 mm - 1 pc.

27. Side rail (drainage hook holder) - 2 pcs.

28. Drainage hooks - 8 pcs.

29. Built-in shelf for linen (if necessary) - 1 pc.

- 30. Multifunctional table (if necessary) - 1 pc.
- 31. Bottom lighting under the bed (if necessary) - 1 pc.
- 32. Warning signal for activating the movement of the bed (if necessary) - 1 pc.
- 33. Belt holders for patient fixation - 4 pcs.
- 34. Operation manual - 1 pc.

1.1.1. Serial numbers: 263437 - 267436

1.2. Medical functional electric 4-section NewCare S bed consisting of:

- 1. Bed frame - 1 pc.
- 2. Bed platform, composition:
 - removable sections of the stock made of HPL-plastic - 3 pcs.
 - fixed section of the stock made of steel - 1 pc.
- 3. Plastic side rails of the back section (if necessary), execution options:
 - with built-in control panel (if necessary) - 2 pcs.
 - with built-in backlit control panel (if necessary) - 2 pcs.
 - without control panel (if necessary) - 2 pcs.
- 4. Plastic side rails of the leg section (if necessary) - 2 pcs.
- 5. Metal side rails, execution options:
 - stainless steel (if necessary) - 2 pcs.
 - painted steel (if necessary) - 2 pcs.
- 6. Head end, execution options:
 - removable with hand grips (if necessary) - 1 pc.
 - removable straight (if necessary) - 1 pc.
 - non-removable (if necessary) - 1 pc.
- 7. Foot end, execution options:
 - removable with hand grips with built-in control panel (if necessary) - 1 pc.
 - removable with hand grips without a control panel (if necessary) - 1 pc.
 - removable straight line with built-in control panel (if necessary) - 1 pc.
 - removable straight line without control panel (if necessary) - 1 pc.
- 8. Wheels, execution options:
 - metal 150 mm (if necessary) - 4 pcs.
 - integral 150 mm (if necessary) - 4 pcs.
 - double 150 mm (if necessary) - 4 pcs.
- 9. Central brake, execution options:

- multifunctional pedal (if necessary) - 1 pc.
 - side pedals of the central brake on the wheels at the foot end (if necessary) - 2 pcs.
 - side pedals of the central brake for all wheels (if necessary) - 4 pcs.
10. 5th guide wheel 100 mm (if necessary) - 1 pc.
11. Foot control pedal (if necessary) - 1 pc
12. Anti-decubitus mattress, variants:
- Density 100 (if necessary) - 1 pc.
 - Density 150 (if necessary) - 1 pc.
 - Density 250 fresh (if necessary) - 1 pc.
 - Density 400 (if necessary) - 1 pc.
 - Density 450 (if necessary) - 1 pc.
 - Density 475 fresh (if necessary) - 1 pc.
 - Density 600 (if necessary) - 1 pc.
 - Density 650 fresh (if necessary) - 1 pc.
 - Density 900 (if necessary) - 1 pc.
13. Hygienic mattress cover, variants:
- Cugar (if necessary) - 1 pc.
 - Select (if necessary) - 1 pc.
14. Additional mattress for extending the bed with a hygienic cover (if necessary) - 1 pc.
15. Lever for emergency lowering of the back section - 2 pcs.
16. Protective bumpers, execution options:
- diameter 84 mm (if necessary) - 4 pcs.
 - diameter 104 mm (if necessary) - 4 pcs.
17. X-ray cassette holder with a guide (if necessary) - 1 pc.
18. Electric drive - 4 pcs.
19. Pneumatic leg section adjustment - 1 pc.
20. Mattress holders - not more than 4 pcs.
21. Mounting holes - not more than 4 pcs.
22. Control panel for medical personnel (if necessary) - 1 pc.
23. Patient control panel, execution options:
- backlit, without blocking of functions (if necessary) - 1 pc.
 - without blocking of functions (if necessary) - 1 pc.
 - with blocking of functions (if necessary) - 1 pc.
24. Power cable - 1 pc.
25. Battery (if necessary) - 1 pc.

26. Built-in segment for lengthening the bed 300 mm - 1 pc.
27. Side rail (drainage hook holder) - 2 pcs.
28. Drainage hooks - 8 pcs.
29. Built-in shelf for linen (if necessary) - 1 pc.
30. Multifunctional table (if necessary) - 1 pc.
31. Bottom lighting under the bed (if necessary) - 1 pc.
32. Warning signal for activating the movement of the bed (if necessary) - 1 pc.
33. Belt holders for patient fixation - 4 pcs.
34. Operation manual - 1 pc.

1.2.1. Serial numbers: 267437- 268436

1.3. Medical functional electric 4-section NewCare XS bed consisting of:

1. Bed frame - 1 pc.
2. Bed platform, composition:
 - removable sections of the stock made of HPL-plastic - 3 pcs.
 - fixed section of the stock made of steel - 1 pc.
3. Plastic side rails of the back section (if necessary), execution options:
 - with built-in control panel (if necessary) - 2 pcs.
 - with built-in backlit control panel (if necessary) - 2 pcs.
 - without control panel (if necessary) - 2 pcs.
4. Plastic side rails of the leg section (if necessary) - 2 pcs.
5. Metal side rails, execution options:
 - stainless steel (if necessary) - 2 pcs.
 - painted steel (if necessary) - 2 pcs.
6. Head end, execution options:
 - removable with hand grips (if necessary) - 1pc.
 - removable straight (if necessary) - 1pc.
 - non-removable (if necessary) - 1pc.
7. Foot end, execution options:
 - removable with hand grips with built-in control panel (if necessary) - 1 pc.
 - removable with hand grips without a control panel (if necessary) - 1 pc.
 - removable straight line with built-in control panel (if necessary) - 1pc.
 - removable straight line without control panel (if necessary) - 1pc.
8. Wheels, execution options:
 - metal 150 mm (if necessary) - 4 pcs.

- integral 150 mm (if necessary) - 4 pcs.
 - double 150 mm (if necessary) - 4 pcs.
9. Central brake, execution options:
- multifunctional pedal (if necessary) - 1 pc.
 - side pedals of the central brake on the wheels at the foot end (if necessary) - 2 pcs.
 - side pedals of the central brake for all wheels (if necessary) - 4 pcs.
10. 5th guide wheel 100 mm (if necessary) - 1 pc.
11. Foot control pedal (if necessary) - 1 pc
12. Anti-decubitus mattress, variants:
- Density 100 (if necessary) - 1 pc.
 - Density 150 (if necessary) - 1 pc.
 - Density 250 fresh (if necessary) - 1 pc.
 - Density 400 (if necessary) - 1 pc.
 - Density 450 (if necessary) - 1 pc.
 - Density 475 fresh (if necessary) - 1 pc.
 - Density 600 (if necessary) - 1 pc.
 - Density 650 fresh (if necessary) - 1 pc.
 - Density 900 (if necessary) - 1 pc.
13. Hygienic mattress cover, variants:
- Cugar (if necessary) - 1 pc.
 - Select (if necessary) - 1 pc.
14. Additional mattress for extending the bed with a hygienic cover (if necessary) - 1 pc.
15. Lever for emergency lowering of the back section - 2 pcs.
16. Protective bumpers, execution options:
- diameter 84 mm (if necessary) - 4 pcs.
 - diameter 104 mm (if necessary) - 4 pcs.
17. X-ray cassette holder with a guide (if necessary) - 1 pc.
18. Electric drive - 4 pcs.
19. Pneumatic leg section adjustment - 1 pc.
20. Mattress holders - not more than 4 pcs.
21. Mounting holes - not more than 4 pcs.
22. Control panel for medical personnel (if necessary) - 1 pc.
23. Patient control panel, execution options:
- backlit, without blocking of functions (if necessary) - 1 pc.
 - without blocking of functions (if necessary) - 1 pc.

- with blocking of functions (if necessary) - 1 pc.
- 24. Power cable - 1 pc.
- 25. Battery (if necessary) - 1pc.
- 26. Built-in segment for lengthening the bed 300 mm - 1 pc.
- 27. Side rail (drainage hook holder) - 2 pcs.
- 28. Drainage hooks - 8 pcs.
- 29. Built-in shelf for linen (if necessary) - 1 pc.
- 30. Multifunctional table (if necessary) - 1 pc.
- 31. Bottom lighting under the bed (if necessary) - 1 pc.
- 32. Warning signal for activating the movement of the bed (if necessary) - 1 pc.
- 33. Belt holders for patient fixation - 4 pcs.
- 34. Operation manual - 1 pc.

1.3.1. Serial numbers: 268437 - 269436

1.4. Medical functional electric 4-section NewCare XSS bed consisting of:

1. Bed frame - 1 pc.
2. Bed platform, composition:
 - removable sections of the stock made of HPL-plastic - 3 pcs.
 - fixed section of the stock made of steel - 1 pc.
3. Plastic side rails of the back section (if necessary), execution options:
 - with built-in control panel (if necessary) - 2 pcs.
 - with built-in backlit control panel (if necessary) - 2 pcs.
 - without control panel (if necessary) - 2 pcs.
4. Plastic side rails of the leg section (if necessary) - 2 pcs.
5. Metal side rails, execution options:
 - stainless steel (if necessary) - 2 pcs.
 - painted steel (if necessary) - 2 pcs.
6. Head end, execution options:
 - removable with hand grips (if necessary) - 1pc.
 - removable straight (if necessary) - 1pc.
 - non-removable (if necessary) - 1pc.
7. Foot end, execution options:
 - removable with hand grips with built-in control panel (if necessary) - 1 pc.
 - removable with hand grips without a control panel (if necessary) - 1 pc.
 - removable straight line with built-in control panel (if necessary) - 1pc.

- removable straight line without control panel (if necessary) - 1 pc.
8. Wheels, execution options:
- metal 150 mm (if necessary) - 4 pcs.
 - integral 150 mm (if necessary) - 4 pcs.
 - double 150 mm (if necessary) - 4 pcs.
9. Central brake, execution options:
- multifunctional pedal (if necessary) - 1 pc.
 - side pedals of the central brake on the wheels at the foot end (if necessary) - 2 pcs.
 - side pedals of the central brake for all wheels (if necessary) - 4 pcs.
10. 5th guide wheel 100 mm (if necessary) - 1 pc.
11. Foot control pedal (if necessary) - 1 pc
12. Anti-decubitus mattress, variants:
- Density 100 (if necessary) - 1 pc.
 - Density 150 (if necessary) - 1 pc.
 - Density 250 fresh (if necessary) - 1 pc.
 - Density 400 (if necessary) - 1 pc.
 - Density 450 (if necessary) - 1 pc.
 - Density 475 fresh (if necessary) - 1 pc.
 - Density 600 (if necessary) - 1 pc.
 - Density 650 fresh (if necessary) - 1 pc.
 - Density 900 (if necessary) - 1 pc.
13. Hygienic mattress cover, variants:
- Cugar (if necessary) - 1 pc.
 - Select (if necessary) - 1 pc.
14. Additional mattress for extending the bed with a hygienic cover (if necessary) - 1 pc.
15. Lever for emergency lowering of the back section - 2 pcs.
16. Protective bumpers, execution options:
- diameter 84 mm (if necessary) - 4 pcs.
 - diameter 104 mm (if necessary) - 4 pcs.
17. X-ray cassette holder with a guide (if necessary) - 1 pc.
18. Electric drive - 4 pcs.
19. Pneumatic leg section adjustment - 1 pc.
20. Mattress holders - not more than 4 pcs.
21. Mounting holes - not more than 4 pcs.
22. Control panel for medical personnel (if necessary) - 1 pc.

23. Patient control panel, execution options:

- backlit, without blocking of functions (if necessary) - 1 pc.
- without blocking of functions (if necessary) - 1 pc.
- with blocking of functions (if necessary) - 1 pc.

24. Power cable - 1 pc.

25. Battery (if necessary) - 1pc.

26. Built-in segment for lengthening the bed 300 mm - 1 pc.

27. Side rail (drainage hook holder) - 2 pcs.

28. Drainage hooks - 8 pcs.

29. Built-in shelf for linen (if necessary) - 1 pc.

30. Multifunctional table (if necessary) - 1 pc.

31. Bottom lighting under the bed (if necessary) - 1 pc.

32. Warning signal for activating the movement of the bed (if necessary) - 1 pc.

33. Belt holders for patient fixation - 4 pcs.

34. Operation manual - 1 pc.

1.4.1. Serial numbers: 269437 - 270436

Accessories*:

1. II-pole (Infusion stand with two hooks) - not more than 4 pcs.
2. IV-pole (Infusion stand with four hooks) - not more than 4 pcs.
3. Oxygen cylinder holder - not more than 4 pcs.
4. Trapeze bar, non-adjustable - not more than 4 pcs.
5. Trapeze bar, adjustable - not more than 4 pcs.
6. Traction frame (Balkan frame) - not more than 1 piece.

*** Accessories are common to all bed variants.**

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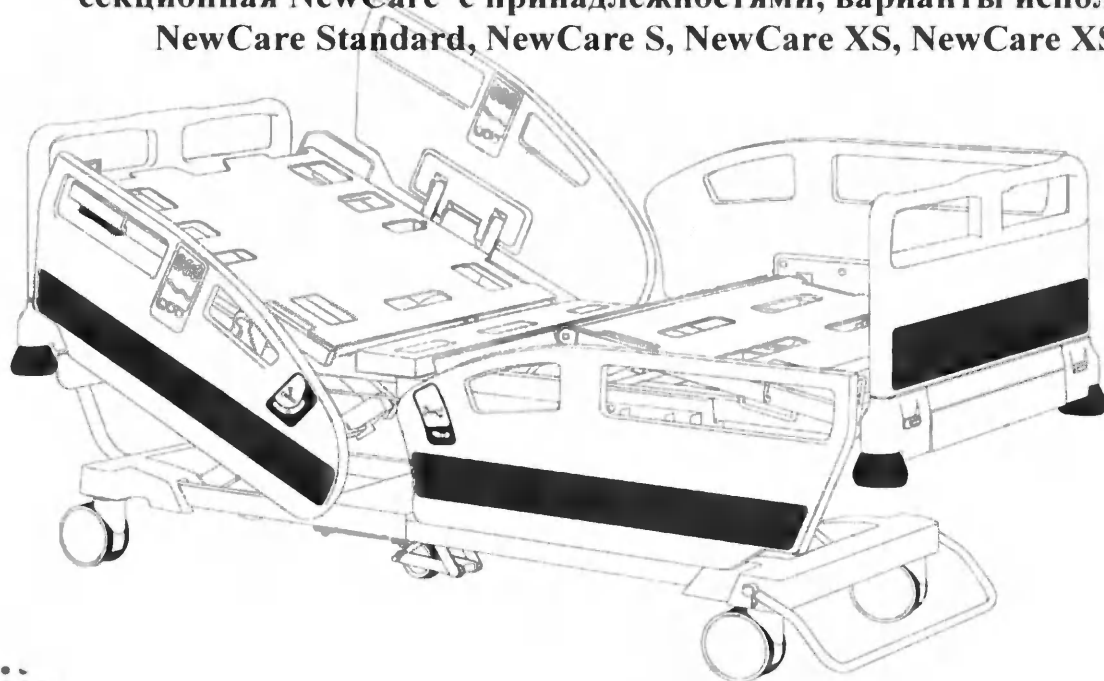
- NewCare Standard
- NewCare S
- NewCare XS
- NewCare XSS

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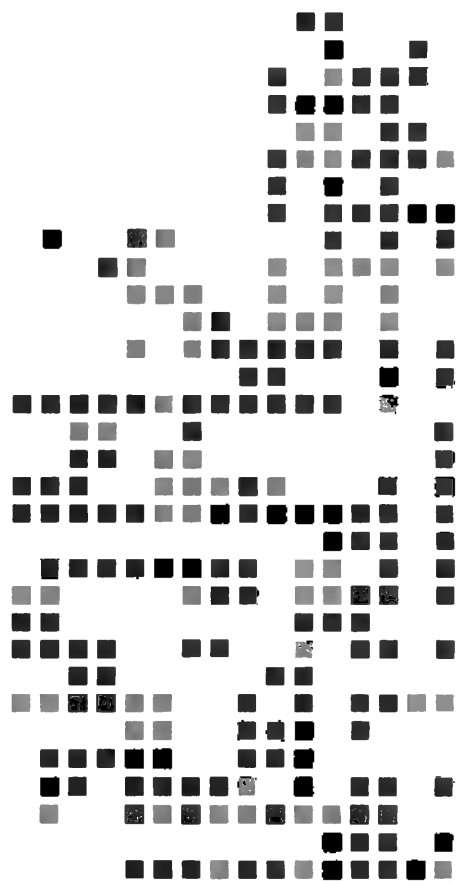
Tel: +34 976 300 033 - Fax: +34 976 320 647

Index

SYMBOLS, LABELS AND MANUAL ORGANIZATION	6
Typographical conventions	7
Symbology	7
Organization	7
Marking	8
SAFETY ADVICES	9
General notices	10
Electrical risk warnings	11
Entrapment warnings	11
REGULATIONS, RULES AND CODES	12
DESCRIPTION	14
Bed Description	15
Bed components	17
Electric components	17
Controls	17
Control box	17
Battery	17
Linear electric actuators	18
Mechanical components	19
Rails	19
Headboard and footboard	19
Weels	19
Centralised braking bar	20
Feet module universal positioning	20
CPR bilateral actioning system	20
Other components	21
Guides and X-ray cassette carrier	21
Protection discs	21
Accessories support	21
Multifunction table	21
Accessories	21
Infusion stands	21
Trapeze	21
Aluminium Balkan frame	21
Multifunction table	21
Tank holder	21
Turn on / off the light under the bed	20
Warning signal	20
Components	22
BED OPERATION	24
Reception	25
Installation	25
Initial operation	25
Controls functioning	26
Control types	26
Wired patient control	26
Rails membrane control	28
Nurse controls	28
Operating	30
Basic control operating	30
Pre-programmed functions	31

Light indicators	32
Lock of function	32
Other functions	33
Bed components functioning	34
CPR bilateral actioning system	34
Rails	34
Plastic full length rails	34
Metal rails	35
Removable Headboard / Footboard with locking system	35
Placement	35
Removal	36
Multifunction pedal	37
Bed braking	37
Bed moving	37
Directional wheel actioning	37
Fifth wheel	38
Feet module positioning system	38
Cassette carrier	39
Bed base panels	40
Movement indicator	40
CLEANING AND MAINTENANCE	41
Bed cleaning	42
Bed maintenance	43
Preventive maintenance guide	44
Quality assurance	44
Product disposal	43
TECHNICAL SPECIFICATIONS	49
Physical bed description	51
Rails physical description	51
Recommended mattress description	52
Use environment (according to direction EN 60601-2-52:2010)	50
Use environment physical description	50
Storage environment physical description	50
Electric system description	52
Electrical diagram	53
“EC” DECLARATION OF CONFORMITY	55
PREVENTIVE MAINTENANCE CONTROL LIST	57
Application	
Mattresses. Technical specification	
Models	
Density 100	
Density 150	
Density 250 fresh	
Density 400	
Density 450	
Density 475 fresh	
Density 600	
Density 650 fresh	
Density 900	

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- **General description.**
Provides a general description of the product to know its main aspects.
- **Functioning and handling**
Describes the use of different mechanical and electrical systems mounted in this product.
- **General procedures, cleaning and maintenance.**
Describes the cleaning tasks and preventive maintenance that have to be taken for a correct product running.

Marking on the device

	Read the user manual
Input: (120V-220V) - 50/60Hz 700W	Motor voltage
	Electric insulation Class II
	Protection against electric hazard accidents - Type B
	IP66 Protection
	Suitable for indoor use only
	Safety Weight Load (SWL)
	Maximum Patient Weight (MPW)
	Do not load over 140 kg. at the beds ends whilst it is extended
	Serial number label and CE mark
	Maximum load weight warning for the under footboard auxiliary table
	Do not store any object under the bed
	Entrapment danger warning in these areas whilst lowering the bed
	Manual emergency CPR
	Recommended mattress dimensions
	Do not throw away the device along with the domestic waste. Recyclable product
	Do not seat on the articulated bed base panels whilst are raised, neither activate movements if there are people seated on the modules

 	Medical device conformity mark to European Union standards
Symbols and labels on the product	
	Catalog number
	Serial number
	Temperature range
 	Do not dispose of equipment with household waste. Disposal must comply with local regulations. Please refer to the Recycling section for details
	Eurasian Conformity

Note: Some symbols are not found on all equipment.

Symbols	Description
	Do not transport on a cart
	Keep away from moisture

	Industrias Hidraulikas Pardo SL Company 's logotype
	Fragile. Handle with care
	Stacking is prohibited

Safety advices

General notices

- This bed is not designed to be used with special anatomy users (amputees, etc) neither with children, because his anatomy or special circumstances can cause entrapment risks not contemplated in our risks tests. This is the reason why this product is not recommended for such cases.
- This manual describes all beds functions as well as its use, explaining all functions and components in detail. To obtain more accurate information relative to the bed handling refer to the manuals specific sections of the manual.
- Refer to the manual before using the bed. Otherwise personal injuries or material damages may occur.
- The determination and problem solutions of the bed can be only done by qualified and authorised technical staff. The determination and problem solutions of the bed done by non authorised personnel can lead to personal injuries or material damages.
- Do not remove RFID labels from the bed. The reference that appears on it are essential to identified the bed and warranties the safety maintenance operations. Otherwise personal injuries or material damages may occur.
- Never raise the bed with other tool than the designed by the manufacturer. Not complying with this instruction may cause a rising mechanism dysfunction and body injuries.
- The motorised mechanisms of the machine can cause injuries. Do not use the bed controls until all the people are away from it. A bed articulation function can be interrupted releasing the control, pressing the button of the contrary movement or unplugging the supply connector from mains.
- The raising or lowering of the bed sections can be stopped by releasing the button on the control device, pressing the section movement button in the opposite direction, or by disconnecting the bed from the power source.
- The bed is designed for indoor use. Do not use or storage it on outdoors or humid environments.
- Do not use the equipment if the width, height and patient weight exceed the limits specified in this manual. Patient injuries may occur.
- Do not use in this bed mattress that cannot be adapted to the bed base articulations.
- Only authorised maintenance staff should make the preventive maintenance of the bed. A preventive maintenance mad by unauthorised personnel could produce personal injuries or material damages.
- The motors use cycles (intermittent operation) does not exceed 2 minutes / 20 minutes, otherwise equipment damages can be caused.
- Refer to the beds user manual before using the accessories to avoid any patient or user injuries risk.
- Only the following accessories specifically identified can be used in this bed. The use of any other accessories not included in the list will compromised the bed and patient safety.
- Do not hit the different movement controls. In the same way, never manipulate the inside of the manual or motor controls.
- Do not activate any movement before placing the mattress previously.
- Unplug the bed before any cleaning operation, bed moving or electric equipment repair is taking place. There is electrical shock risk.
- Never manipulate the manual or motor control inside.
- The electrical installation which the bed is connected has to have the regulatory protections, otherwise do not connect the device.
- In case any deformation is observed in any of the bed element, get in contact with the distributor.
- In the modules rising and lowering movements, it has to be observed that no materials or body parts are in-between the bed base and the bed supporting frames. Otherwise there are entrapment risks between the bed base and its supporting points whilst lowering.
- For the patient safety, the bed should be in it lowest position when the patient is not being attended, to avoid the risk of falling damages.

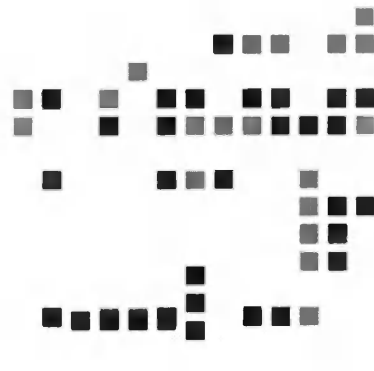
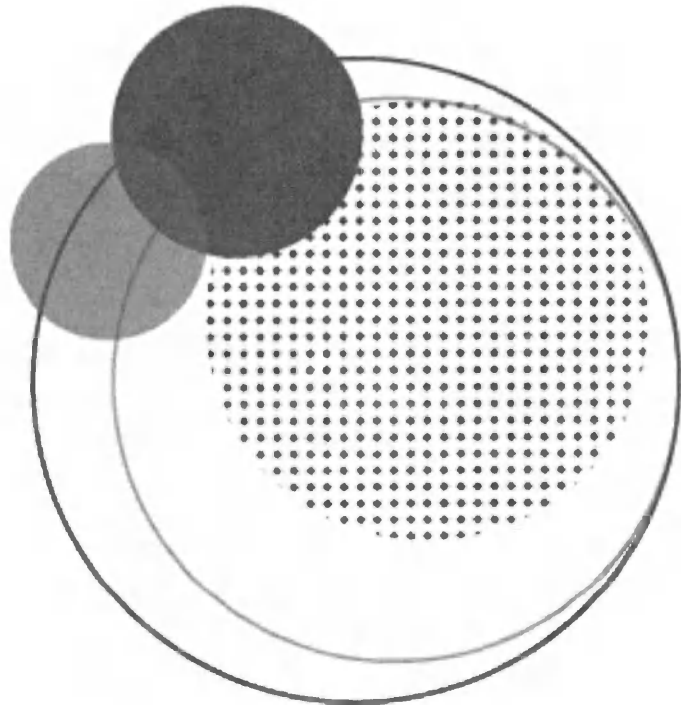
Electrical risk warnings

- The use of electric equipment involves electrical risk. The staff has to be capable and have the respective information about the inherent risks of electrical equipment.
- An abStandardlly leakage of current is symptom of a AC power cable and of power supply unit. If the leakage is over 500 microamperes, personal injuries can be caused.
- Unplug the bed before any bed cleaning or maintenance action is made. Otherwise, personal injuries or material damages may occur. Refer to the specific manual section for further warnings.
- For patient or staff safety, take precaution whilst connecting other wired devices. Avoid harming or pressing the wires of other equipment with the different elements and mobile parts of the bed.
- Concerning the bed connection, avoid harming or pulling the power cord, because it could imply electrical risk.

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Regulations,
rules and
codes

Regulations, rules and codes

NewCare V2 Standard bed is in accordance with the following regulations:

- Directive 93/42/CE
- Directive 2006/42/CE
- Directive 2006/95/CE
- Directive 2004/108/CE
-

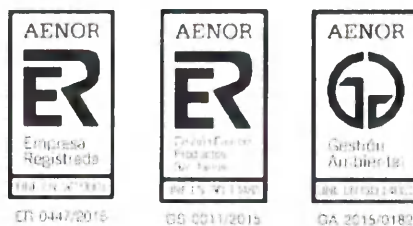
NewCare V2 Standard bed complies with the following rules:

- Directive UNE-EN ISO 12100
- Directive UNE 20324
- Directive UNE-EN 60204-1
- Directive UNE-EN 60601-1
- Directive UNE-EN 60601-1-2
- Directive EN 60601-2-52
-

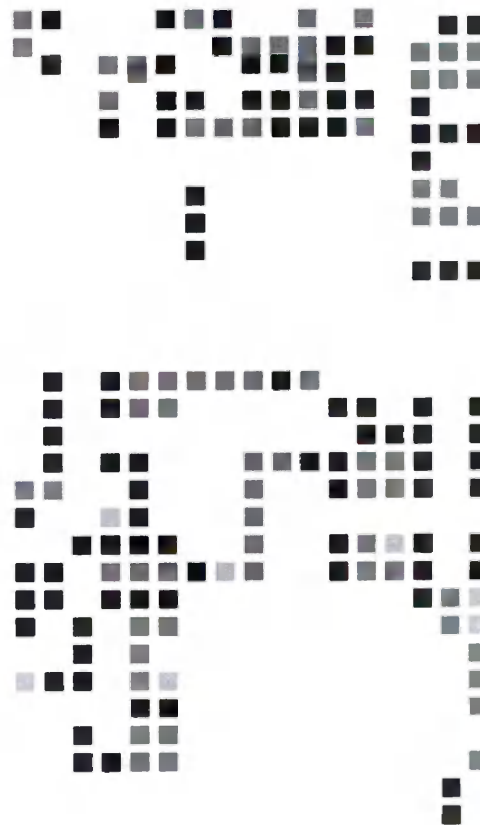
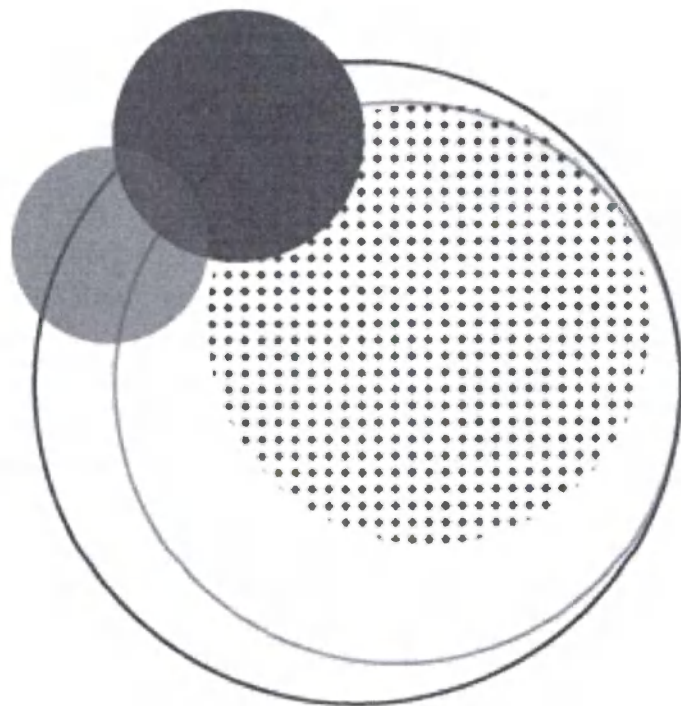
NewCare V2 Standard bed has the following marking:

Industrias Pardo s.l. has the following certifications:

Description

Description

The hospitalization bed NewCare V2 Standard is specially created to be used in hospitals, clinics and other similar uses. Designed for 10 to 15 years lifetime, according to its use.

Its structure is made of steel with texturised coating painting with polyester base. The thickness of the steel is between 1.5 and 3 mm, depending on the area of the product. The coating is resistant to chemical and mechanical stress.

The bed base is divided in two articulated panels electrically actioned (backrest and extremities), another manually articulated panel (feet) and one fix panel. The bed base is made of four HPL (High Pressure Laminated) high resistance pieces, removable and fireproof one fixed panel made of polyester-based fine embossed steel.

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It has four 150mm wheels. the 5th wheel acts as directional. In case there has not 5th wheel, one of the four wheels of the bed will be antistatic and directional.

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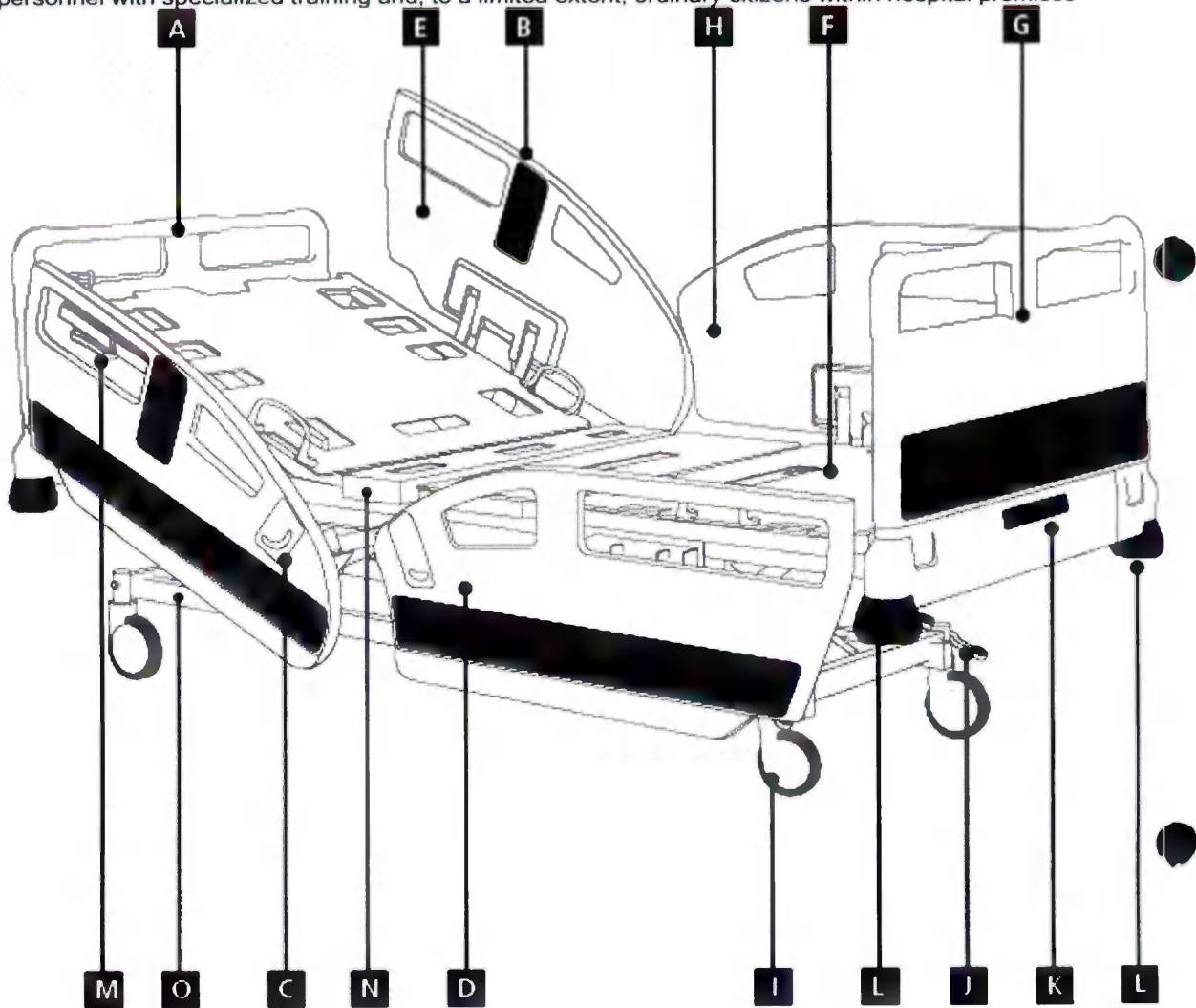
Before using the bed ensure to read and understand the NewCare V2 user manual content. It is important to read and strictly follow the safety instructions contained in the manual.

Intended use

Functional electric 4-section NewCare bed in versions with accessories is designed to care for the patient, increase comfort and create the desired positioning of the patient in a hospital environment.

Potential consumers

- Medical personnel with specialized training and, to a limited extent, ordinary citizens within hospital premises

Элемент	Описание	Элемент	Описание
A	Removable headboard with hand grips	I	Ø 150 mm wheels
B	Integrated control on the side rail	J	Breaking pedal
C	The indicator of the dorsal section inclination	K	The control panel on the footboard
D	Trendelenburg position indicator	L	Protective discs
E	Plastic side rails of the dorsal section	M	Lever for emergency lowering of the dorsal section

F	3 removable bed sections made of HPL	N	Fixed section of the bed made of steel
G	Removable footboard with hand grips	O	Bed frame
H	Plastic side rails of the foot section	P	Mattress Holders

Bed components

Electric components

*Do not hit the controls. Do not manipulate in case of breakage and do not access to the controls inside.
In case of power failure, the control will enter in lock mode. To unlock follow the instructions specified in this manual.
In case of breakdown, please, consult your dealer.*

*The use of electrical equipment carries electric risks. The staff should be capable and must have available the information inherent to the electrical equipment.
The electric circuit is electrically insulated from the bed metal parts.
To avoid motor overload, the movement should not last over two minutes.
Any electric component reparations should be done with the bed unplugged of mains and disconnected battery.
Do not manipulate the boxes, controls or electrical motors in case of breakdown or access inside the electric components.
Observe that the power socket and bed plug are properly connected. An error respecting this can cause overheating in these components.
Check that the mains voltage match with the nominal voltage needed.
The electrical installation where the bed is connected must have the regulatory protections.*

Controls

The bed has different control options depending on the final use and client needs.

NewCare V2 bed is equipped with different controls depending on options and configurations. They have all functions to handle the bed.

- Locking position indicator in nursing control.
- Compatible with right/left-handed.
- Membrane controls integrated at both sides of the side rails.

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Characteristics:

- Supply voltage: 230 V. 50Hz.
 - Operating temperature range: 5° - 40°
- IP66 protection
Power cable

The power cable has a minimum length of 2.5 m from the mains plug. The power cable has a cast plug. Insulated and sheathed power cables and other external flexible bed cables are made of ethylene-propylene rubber and provide adequate protection against damage during normal operation.

Battery

The battery is located in the bed chassis, in the central area next to the extremities motor and control

box.

- Rechargeable battery type: Lead gel.
- Voltage: 24 V
- Capacity: 1,2 Ah.

Cycle: about 300 movements The maximum time a bed with battery can be unplugged is three months, although is recommendable charge it every two months if the bed goes to remain unplugged.

The batteries have a limited life so that it's convenient checking usually its charge condition. It's recommendable replacing them when exhaustion is observed or every three four years. The nurse control has a battery indicator, which shows when the bed is plugged to mains and the battery is properly charging.

Removing the battery: Unscrew bracket ends and removing bolts which are fastening it.

Placing the battery: screw on the bracket ends with bed frame by means of nuts and bolts.

Linear electric actuators

It has 2 linear actuators located in bed chassis which take over with the backrest and extremities sections rising; and 2 actuators in the chart for its lifting and lowering.

- Motor voltage: 24 V
- Installation dimensions:
 - 445 mm. (Backrest).
 - 336 mm. (legs).
 - 363 mm. (bed).
- Thrust:
 - 4000 N. (backrest).
 - 3000 N. (legs).
 - 6000 N. (bed).
- Stroke length:
 - 200 mm. (backrest).
 - 100 mm. (legs).
 - 150 mm. (bed).
- Protection: IP66.

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The side rails are made of high strength polyolefins and are installed on the sides of the bed and attached to the bed base.

Each side rail is equipped with a lever that can be operated with one hand. The lowering of the lever is softened by a gas spring, which prevents it from falling suddenly. The side rails have a safety clip in the upper position and an anti-pinch system that prevents them from locking in the lower position.

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Metal side rails (4 tubes rails in stainless or painted), ordered separately

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Each handrail has a lock and unlock system.

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Raise the leg section only when the patient is in a bent knee position, otherwise there is a risk of injury.

The leg section can be raised or lowered to the desired position using a gas spring.

Changing the position of the leg section of the bed helps patients with medium to long-term stay to change position, and also helps prevent the formation of pressure ulcers in the ankle area.

Headboard and footboard

Removable headboard and footboard with locking system.

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The removable headboard and footboard of the bed are made of polypropylene and have rounded corners and edges, as well as ergonomic holes for a comfortable grip while moving the bed.

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Wheels

NewCare V2 Standard beds have four integral wheels (diameter 150 mm), with high resistance to rolling, wear and rust, as well as with minimal noise during movement.

Other wheel models can be purchased separately: metal wheels with a diameter of 150 mm or 125 mm metal wheels with an individual brake.

Centralised braking bar

*Remember always to keep the brake on when the bed is not being moved, to avoid injuries in case of accidental displacements.
Ensure, at any bed transportation, that it is unplugged from mains.*

Centralised braking side pedals (sold separately)

Placed at the feet side of the bed, the multifunction pedal allows the centralised braking with auto locking of the four wheels, total free movement and a third possibility, depending on the bed configuration the 5h wheel can be activated or the directional wheel instead, once they are activated the bed's manoeuvrability will be improved whilst it is being moved, avoiding lateral displacements.

Individual brake pedals (sold separately)

NewCare V2 Standard beds can be equipped with individual brake pedals by selecting wheels compatible with this system. Pedals are designed for centralized braking with automatic locking of all four wheels, for centralized unlocking.

It is possible to purchase an additional warning signal about the movement activation.

Back section emergency lowering lever

Never use the emergency lowering levers of the back section to raise the back section of the bed, as this could damage the lift motor and damage the bed.

The emergency lowering levers of the back section located on both sides of the bed are designed to quickly lower the back section of the bed to a horizontal position.

Holders for accessories

Holders are intended for use with Industrias Pardo accessories only. Use of third-party accessories with holders may result in equipment damage and patient injury.

Around the perimeter of the bed, there are various holders for accessories, such as:

- Mounting holes for infusion stands for intravenous infusion, devices for pulling the patient and Balkan frame;
- Side rail (drain hook holder) located under the bed with adjustable hooks;
 - Belt holders for restraining the patient, located along the sections of the bed.

X-ray cassette guide

The bed has guides in the backrest section to fit the X-ray cassette carrier to be able to make radiographies on the thorax section to the patient without moving him from the bed. There more, the bed base material is radio transparent.

Protection discs

Check the correct installation of the accessories regularly.

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Model	Protection disc	Diameter	Image
Standard	Blue (wide in scope)	Ø104 mm	
S	White (narrow in scope)	Ø84 mm	
XS	Blue (wide in scope)	Ø104 mm	
XSS	White (narrow in scope)	Ø84 mm	

Multifunction table

The multifunctional table is installed on the footboard of the bed using special hooks without the use of additional tools; the hooks provide high reliability and stability of the structure.

The table can be used by both the service staff and the patient independently, for example, during meals or for other purposes.

The maximum load on the table is 5 kg.

Guides and X-ray cassette carrier

The guides installed in the dorsal section of the bed allow you to install a holder for an X-ray cassette. The holder is made of stainless steel and equipped with a handle for which it is convenient to extend it; the height of the holder is adjustable.

The holder for the X-ray cassette can be removed from the dorsal section of the bed without removing the section itself.

Lower illumination under the bed

The bed is equipped with a built-in lighting under the frame, which makes it easier for the attendant to move the bed through dark rooms and is convenient for the patient to use at night to get up and down on the bed safely. The backlight can be switched on at any time using the built-in control panel on the side rail, as well as using the patient control panel with backlight activation function.

Bed motion activation warning signal

NewCare beds are equipped with a bed motion activation warning signal. The signal is triggered when the bed is connected to the power supply and the lock is released from the brake system.

Accessories (sold separately)

Infusion stands

The infusion stand fits into one of the insertion holes; the column consists of two pieces of stainless steel or coated steel, a handle for adjusting the height and an anti-leakage system on the down pipe of the column.

Stands with two or four hooks for infusion bags are available, maximum load per hook is 2 kg.

Unregulated trapeze bar

The trapeze can be hold in one of the supports and it is composed of one steel arm with ergonomic designed handle and a belt with anti-sliding system. The belt length is 300 mm.

Height adjustable trapeze bar

The trapeze can be hold in one of the supports and it is composed of one steel arm with ergonomic designed handle and a belt with anti-sliding system. The belt length is 300 mm and can be adjusted by pressing button.

Aluminium Balkan frame

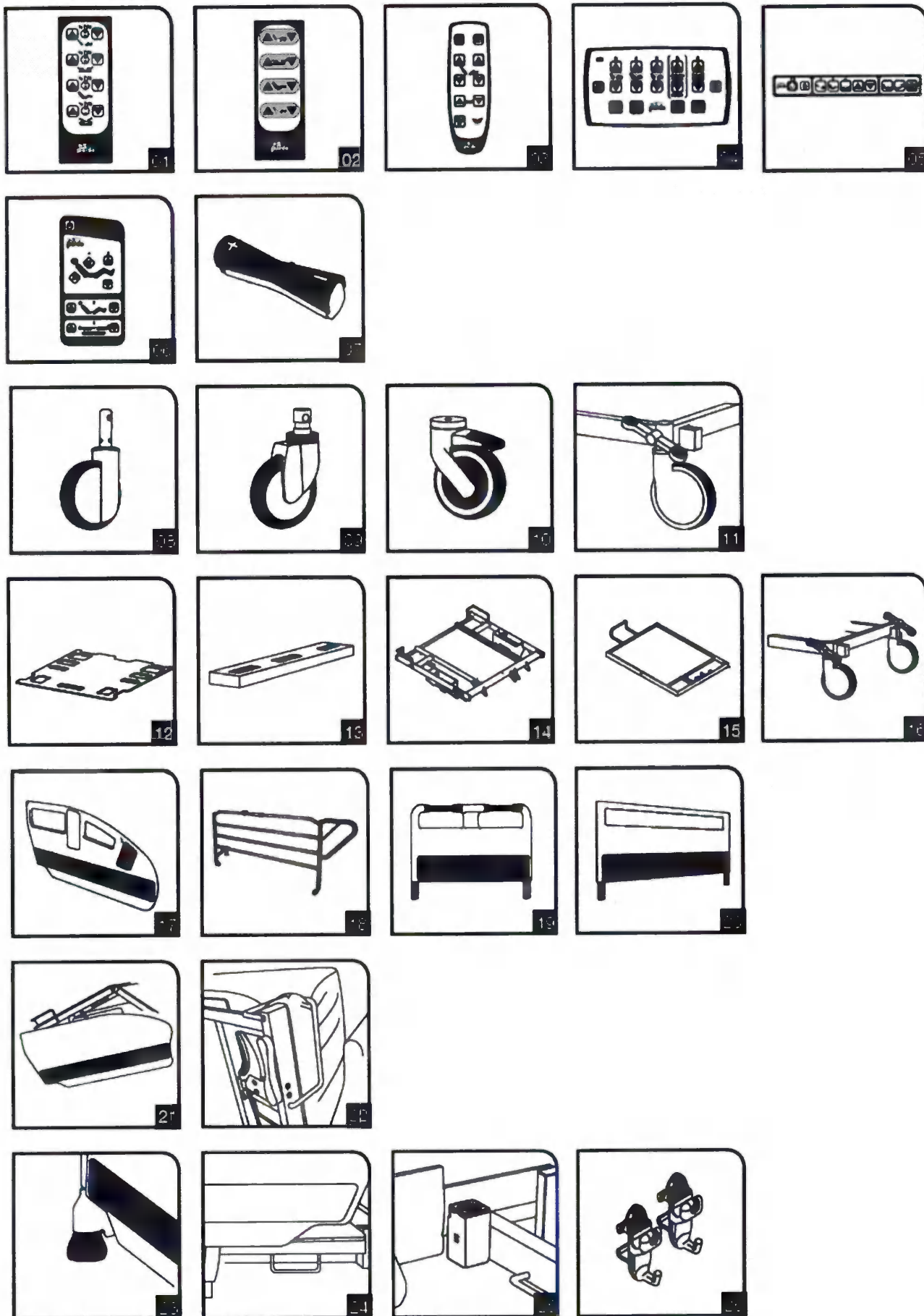
It is a frame structure made of aluminum with an octagonal section.

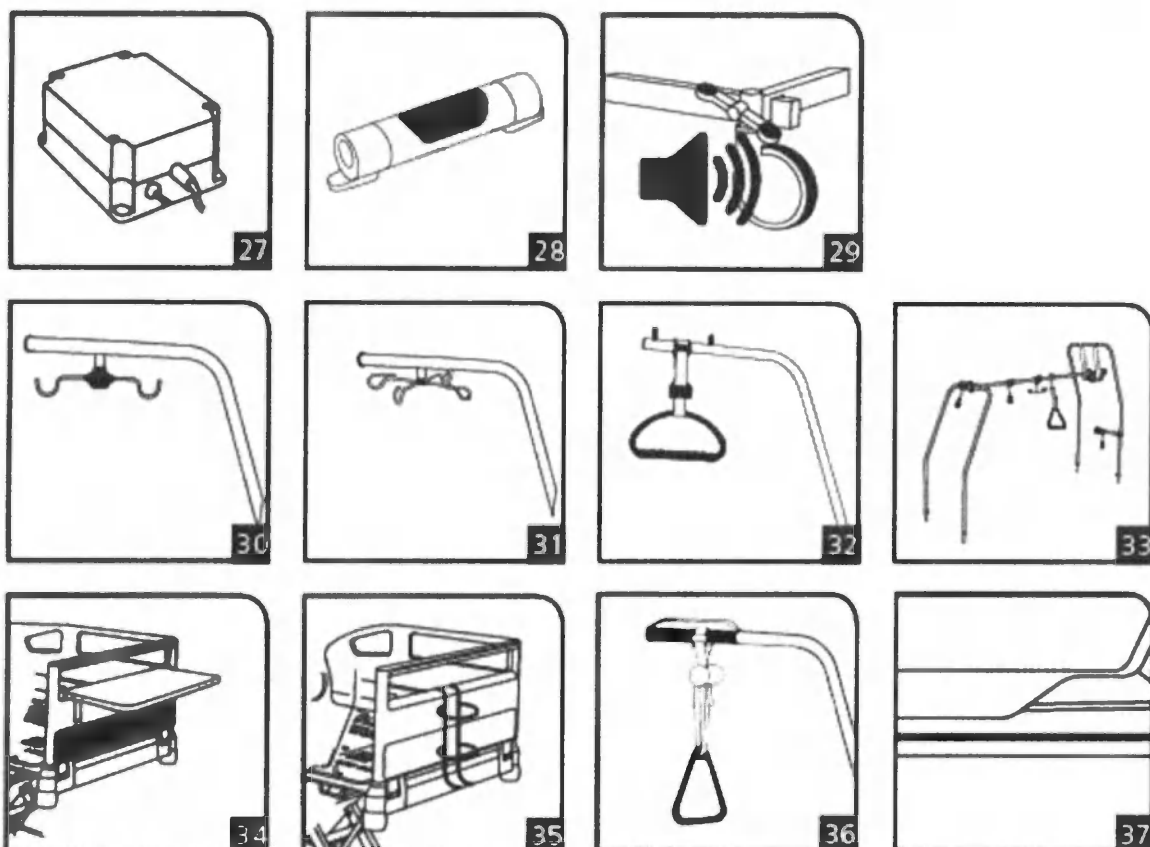
Balkan frame consists of two curved bows secured on both sides of the headboard and footboard of the bed, and one center bar along the bed, fitted with holders, balancers, pull-up device (with handle and strap) and infusion stand. The safe working load is 75 kg.

Tank holder

Metallic structure where can be attached to the footboard or the headboard of the bed and can hold an oxygen bottle. The safe working load is 15 kg.

Components

01	Patient control with functional lock	14	X-ray cassette guide	26	Drainage hooks
02	Basic patient control without functional lock	15	X-ray cassette carrier	27	Battery
03	Light patient control without functional lock	16	Two side center braking pedals	28	Lower illumination under the bed
04	Medical personnel control panel	17	Plastic side rail of the backrest section	29	Bed motion activation warning signal
05	Medical personnel control panel on footboard	18	Metal side rail	30	Infusion stands with two hooks
06	Integrated control panel on the plastic side rail of the backrest section	19	Removable footboard/headboard with hand grips	31	Infusion stands with four hooks
07	Foot control pedal	20	Straight removable footboard/headboard	32	Height adjustable trapeze bar
08	Integral wheel, Ø 150 mm	21	Plastic side rail of foot section	33	Balkan frame
09	Metal wheel, Ø 150 mm	22	Backrest section emergency lowering lever	34	Multifunction table
10	Metal wheel, Ø 125 mm with individual braking	23	Protection discs	35	Tank holder
11	Center braking side pedal	24	Patient strap holder	36	Unregulated trapeze bar
12	Removable bed sections made of HPL plastic	25	Installation holes	37	Side rail (drain hook holder)
13	Fixed steel section				

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Installation

- Ensure not to step on any wire whilst moving the bed.
- Ensure the power supply cable is attached to the headboard hook.
- Ensure the wheels are unlocked before moving the bed.
- Move the bed on appropriate surfaces.
- Mount the headboard and footboard on the bed.
- Ensure the power supply cable is not trapped or overstretched when the bed is placed.
- Ensure the plug is correctly connected.

Initial operation

- Our beds are prepared to operate from the outset. The battery is charged, in case, at the time the bed is plugged, the battery is discharged, it must be recharged before putting it into service.
- An initial bed movement testing can be made:
 - Lift the bed to its highest position.
 - Take down the bed to its lowest position.
 - Articulate the bed base sections.
- Check the rails cushioned lowering.
- Check that the headboard and footboard can be released correctly.
- Check the correct function of all the bed controls

*Only authorized maintenance staff of the establishment should do the start up of the beds.
The manufacturer will not be responsible of any bed or battery damage caused because by an inappropriate handling of the bed or components.*

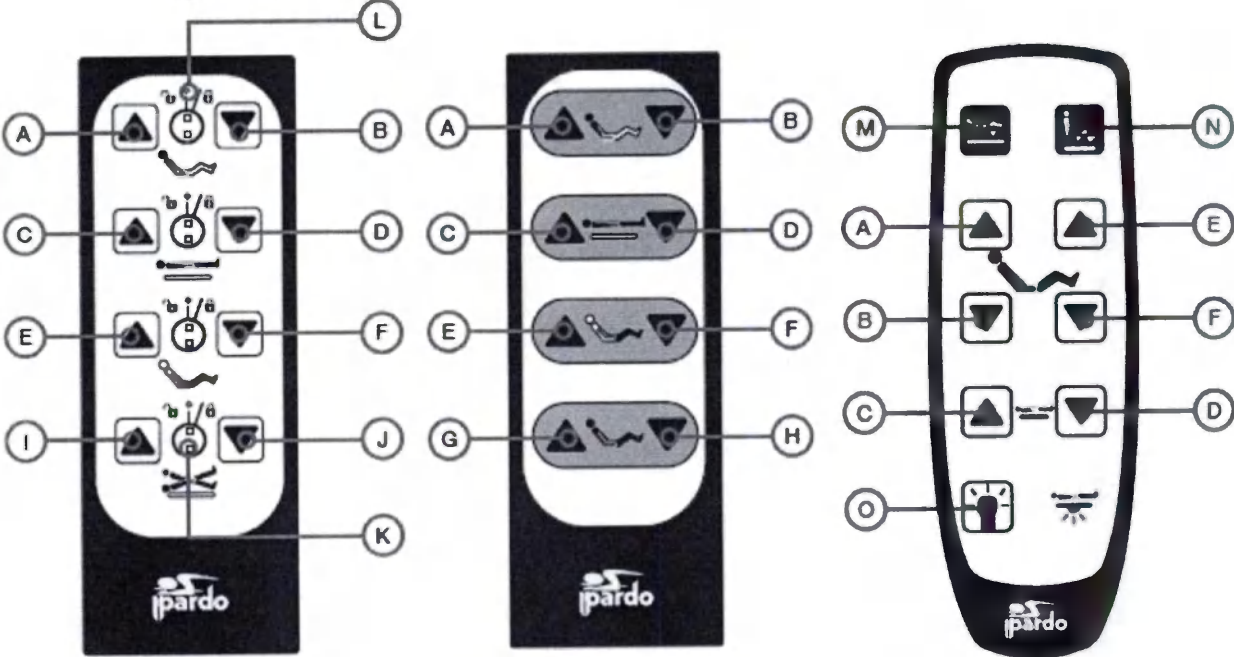
Controls functioning

The bed runs by pressing different control elements. The control elements will depend on the chosen bed model. The set of control panels depends on the bed configuration chosen by the client. The blocking of individual functions, which is only available on the control station for medical personnel, also blocks the selected functions on all other control stations.

Types of control panels

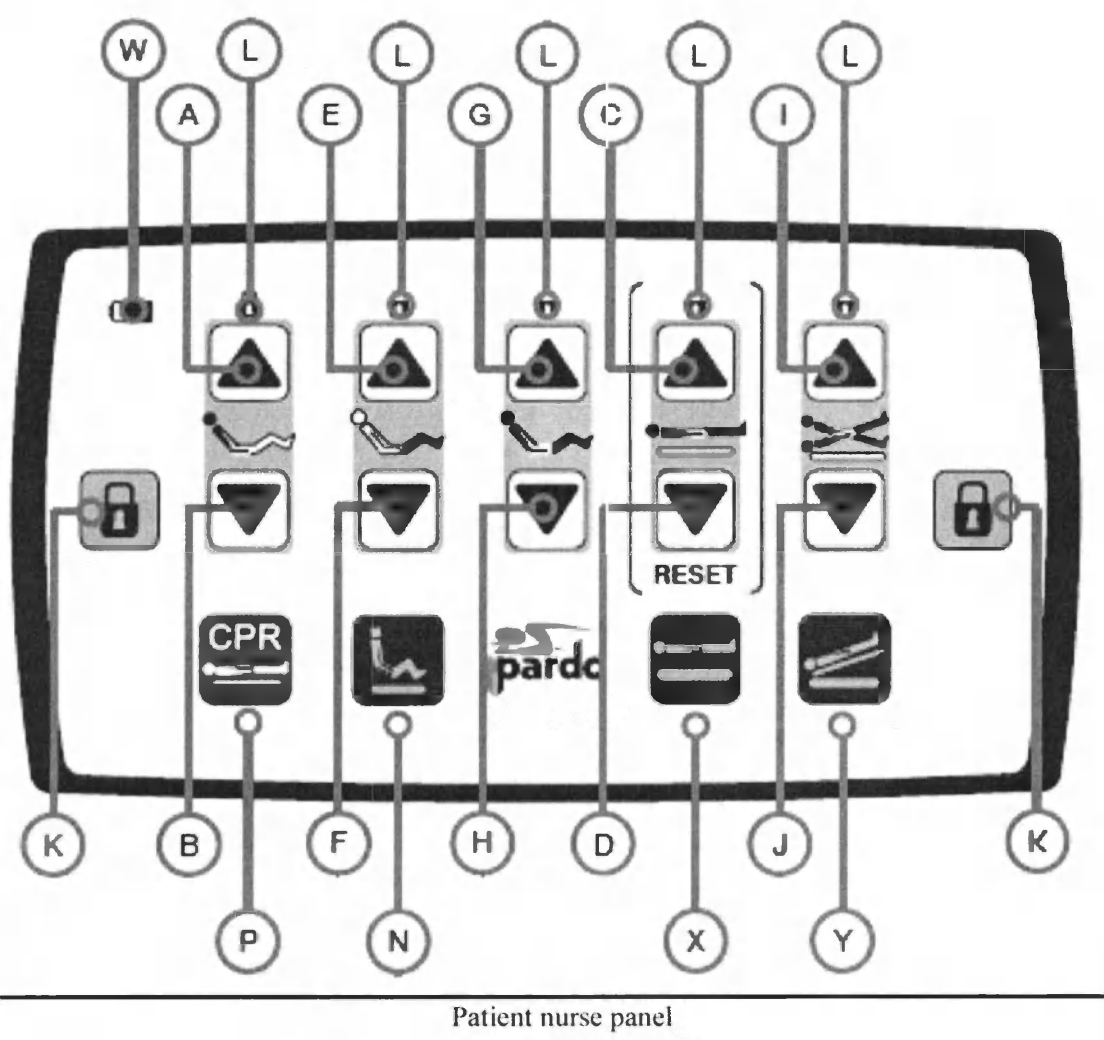
Patient control panel

- Wired patient control:
 - Lock patient control.
 - Basic patient control.
 - Light patient control.

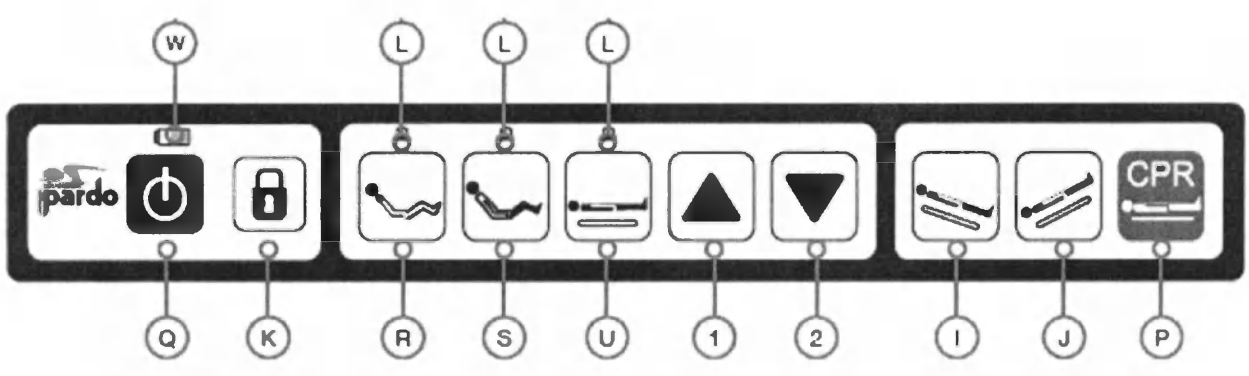
Lock patient control		Basic patient control		Light patient control	
A	Rise backrest section	I	Trendelenburg position		
B	Low backrest section	J	Reverse Trendelenburg position		
C	Rise bed base	K	Lock function		
D	Low bed base	L	Led indicator for locking movement		
E	Rise femoral section	M	Flat		
F	Low femoral section	N	Cardiac chair		
G	Rise backrest and femoral section simultaneously (preset position)	O	Turn on / off the light under the bed		
H	Low backrest and femoral section simultaneously (preset position)				

Patient nurse panel

Patient nurse panel

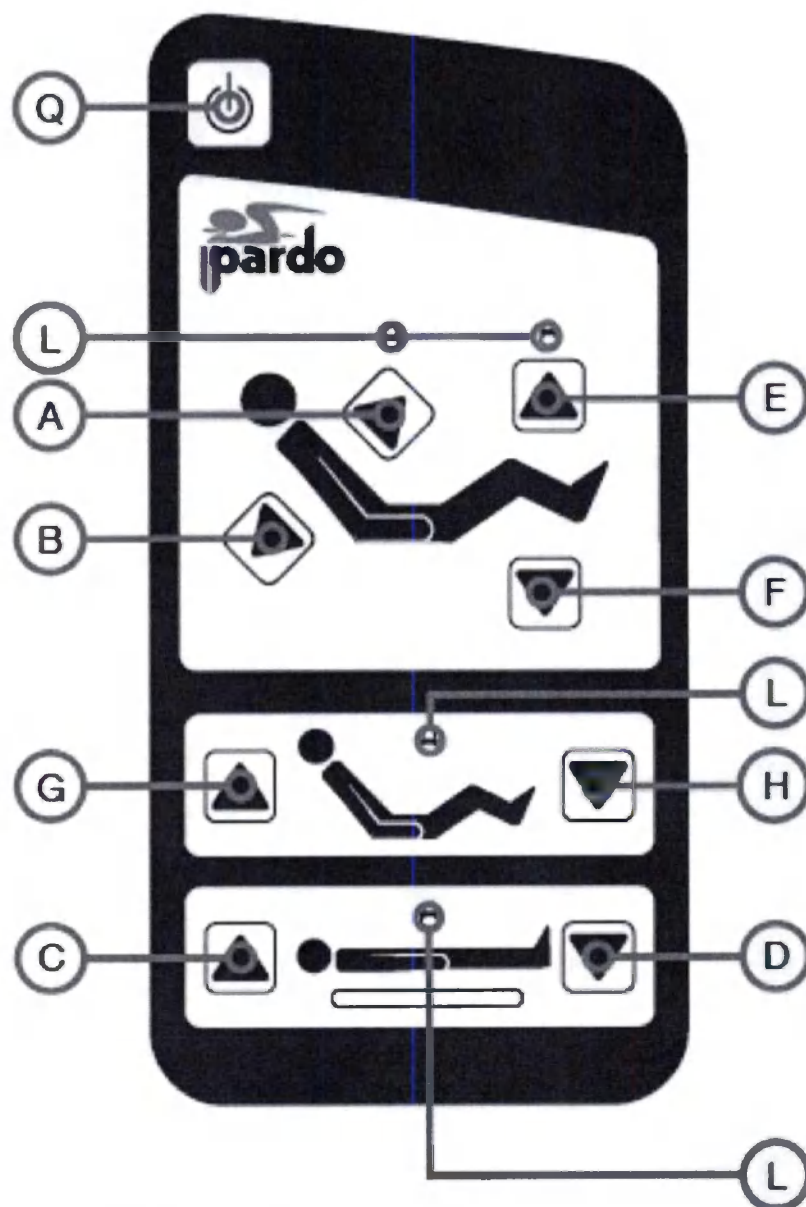
Built-in control panel for nurse on the footboard of the bed

A	Rise backrest section	N	Cardiac Chair
B	Low backrest section	O	Turn on / off the light under the bed
C	Rise bed base	P	Emergency low backrest section
D	Low bed base	Q	ON (activation of control panels) Select backrest

E	Rise femoral section	R	Select backrest section
F	Low femoral section	S	Select femoral section
G	Rise backrest and femoral sections simultaneously (preset position)	U	Select bed base
H	Low backrest and femoral sections simultaneously (preset position)	W	Led indicator for battery charge
I	Trendelenburg position	Y	Emergency Trendelenburg position
J	Reverse Trendelenburg position	X	Inspection position (raises the bed to the highest position)
K	Lock function	1	Rise
L	Led indicator for locking movement	2	Low

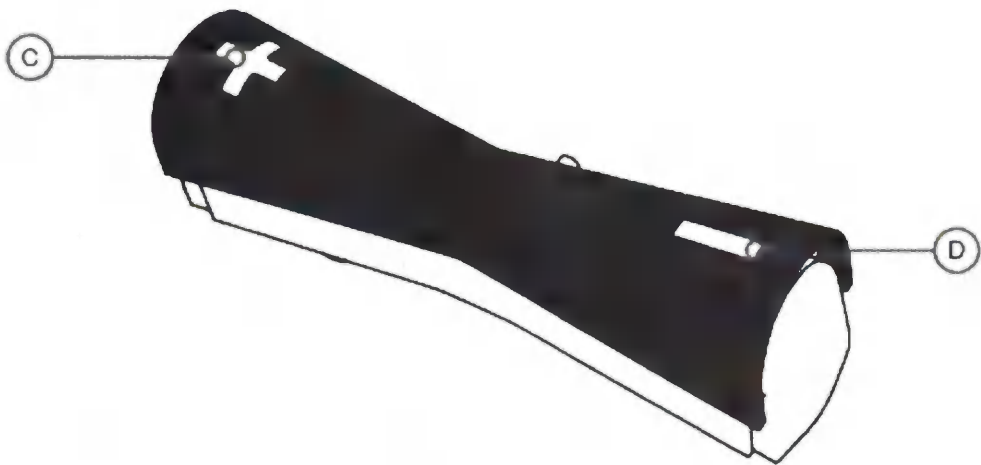
Integrated control panel in the side rail

Integrated control panel in the side rail

A	Rise backrest section	N	Cardiac Chair
B	Low backrest section	F	Low femoral section
C	Rise bed base	G	Rise backrest and femoral sections simultaneously (<i>preset position</i>)
D	Low bed base	H	Low backrest and femoral sections simultaneously (<i>preset position</i>)
E	Rise femoral section	L	Led indicator for locking movement
F	Low femoral section	Q	ON (activation of control panels) Select backrest

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There are crushing risks whilst lowering the bed between it and the cart, in case you are resting your feet on the cart.

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Rise and low backrest section

With the backrest rising and lowering buttons (A and B) the backrest module can be positioned at the desired angle.

To activate the movement, keep on pressing the backrest raising or lowering button to reach the maximum or minimum tilting respectively.

From footboard membrane it's possible doing this action selecting backrest button (R) and pressing raising (1) or lowering arrows (2).

Do not use the RCP bilateral actioning levers to lift the bed's backrest under no circumstances, because it may cause a lifting motor failure and the bed will stop working.

Rise and low femoral section

With the legs rising and lowering buttons (E and F) the legs module can be positioned at the desired angle.

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From footboard membrane it's possible doing this action selecting legs button (S) and pressing raising (1) or lowering arrows (2).

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From footboard membrane it's possible doing this action selecting bed height button (U) and pressing raising (1) or lowering arrows (2).

Raising the bed to the maximum height makes it easier for the staff, while lowering it to the minimum height helps the patient to get up and down on the bed comfortably and safely.

	<i>Do not try to lift or low the bed with other means that the ones designed by the manufacturer. If this warning is not taken in count a lifting mechanism dysfunction or body injuries may occur.</i>
	<i>There are crushing risks whilst lowering the bed between it and the cart, in case you are resting you feet on the cart. There are crushing risks whilst lowering the bed when the rails are all the way down, or if when the rails are being folded with the bed in its lowest position.</i>

Rise or low the backrest and femoral sections simultaneously (preset position)

With the backrest and femoral sections rising and lowering buttons (G and H) the backrest and femoral modules can be positioned at the desired angle at the same time.

To activate the movement keep on pressing the backrest and femoral sections raising or lowering button to reach the maximum or minimum tilting respectively.

Trendelenburg and reverse Trendelenburg

Using the Trendelenburg position (I) and reverse Trendelenburg position buttons (J):

- Trend

Activation: Press the Trendelenburg position button (I), and hold to tilt the foot end of the bed relative to the head end until the maximum tilt angle is reached.

- Reverse Trendelenburg

Activation: Press the reverse Trendelenburg button (J) and hold to tilt the head end of the bed relative to the foot end until the maximum tilt angle is reached. The rest of the bed remains in its original position.

	<i>Trendelenburg and reverse Trendelenburg movements always have to be done under medical supervision.</i>
---	---

Pre-programmed functions

Emergency CPR

Through CPR (P) button, all the bed modules rest in horizontal position and the bed will go down to its lowest position, for a possible urgent cardiopulmonary resuscitation.

	<i>There are entrapments or crushing risks whilst lowering the moving parts of the bed base. Check there are not obstacles on its way down.</i>
---	--

Exploration or maximum height position

Pressing the examination position button (X) raises the bed to the highest position for the comfort of the medical staff.

	<i>There are entrapments or crushing risks whilst lowering the moving parts of the bed base. Check there are not obstacles on its way down.</i>
---	--

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	<i>Trendelenburg and reverse Trendelenburg movements always have to be done under medical supervision.</i>
---	---

Horizontal position

Function available only on the patient control panel. To activate, press the corresponding button (button M) on the control panel, all sections of the bed will lower, and the bed will take a completely horizontal position, convenient for the patient to rest.

Cardiac chair

With the cardiac chair button (N) the modules tilt to its maximum; the bed base tilts to the reverse Trendelenburg position. The movement successions needed to reach this position follows certain order and intervals to improve patient ergonomics and safety.

	<i>There are entrapments or crushing risks whilst lowering the moving parts of the bed base. Check there are not obstacles on its way down.</i>
--	--

On / Off button

By pressing the On / Off (O) button, the rails membrane control is activated to be able to use the functions located in it. After an inactivity of 60 seconds, the membrane turns off and it is necessary to press again the button to reactivate.

	<i>Each membrane has its own activating button. It is necessary to press its own button to activate each membrane.</i>
---	---

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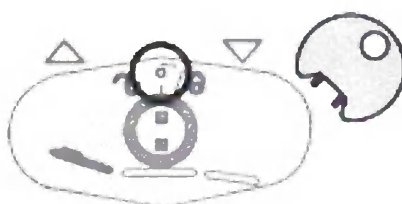
Bed function locking**Patient control with locking**

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Nursing control

Locking different sections or functions will be made pressing the button signalled with a key (K) and then press at the same time the section to be locked (A, C, E, G or I). Indicator light (L) which has been pressed will become fixed. To unlock, redo the same action. The indicator light will remain flickering for a few seconds.

Nursing membranes

- **Backrest movement locking**

Press the indicated locking button with a key (K) and then at the same time press backrest button (R). The indicator light (L) will remain fixed. To unlock, redo the same action. The indicator light will remain flickering for a few seconds.

- **Femoral movement locking**

Press the indicated locking button with a key (K) and then at the same time press backrest button (S). The indicator light (L) will remain fixed. To unlock, redo the same action. The indicator light will remain flickering for a few seconds.

- **Raising/lowering movement locking**

Press the indicated locking button with a key (D) and then at the same time press backrest button (S). The indicator light (L) will remain fixed. To unlock, redo the same action. The indicator light will remain flickering for a few seconds.

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- **Patient control panel**

Keep pressed for some seconds buttons C & D until the beep disappears. The system is reset and restarted to work properly.

- **Control panel for medical personnel**

Keep pressed for some seconds buttons C & D until the beep disappears. The system is reset and restarted to work properly.

- **Control panel for medical personnel at the footboard of the bed**

Keep pressed for some seconds buttons C & D until the beep disappears. The system is reset and restarted to work properly.

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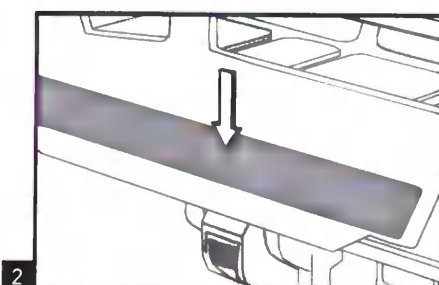
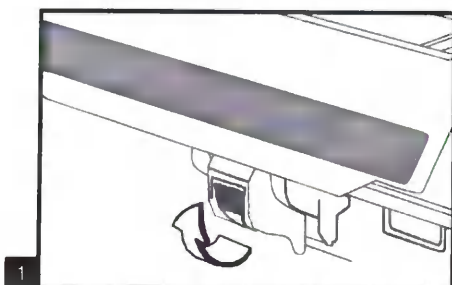
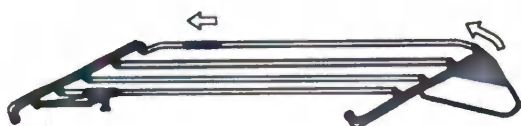
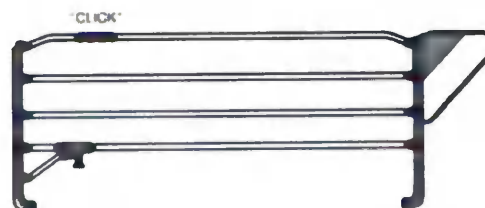
Metal side rails

Side rail raising

- Pull of the side rail (1) towards headboard up to hearing "click" from the lower handle (2), then the side rail is locked in highest position.

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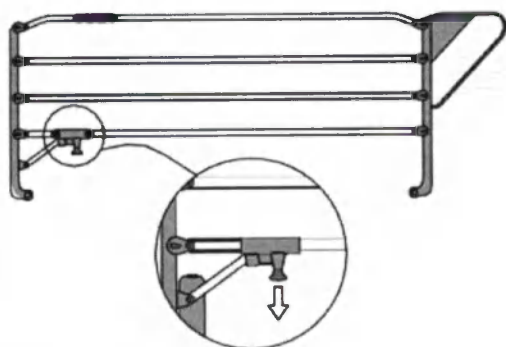
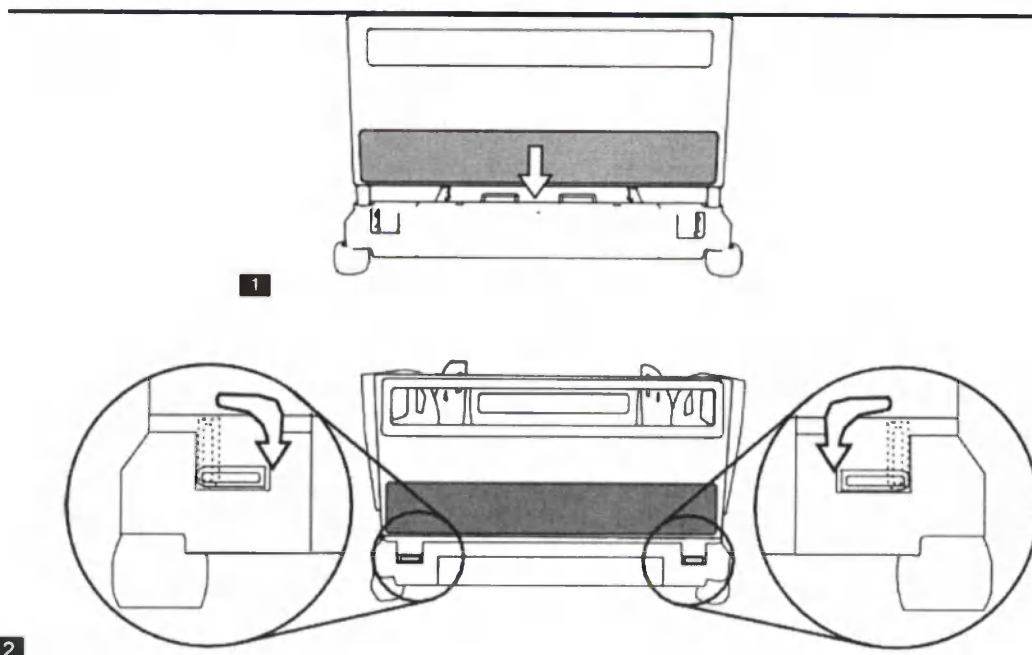
Removable Headboard / Footboard

Placement

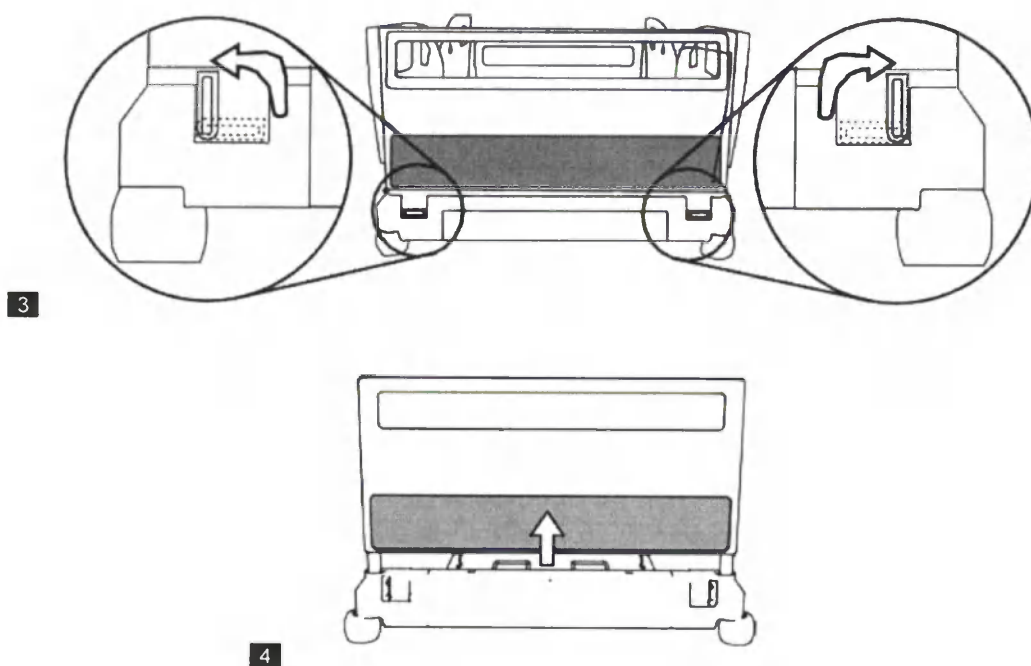
Place the footboard anchors over the beds corner holes. Insert until they are fully in place (1). Lock by turning the levers at each side of the bed towards the inside of the bed until they are horizontal, as shown in the drawing (2).

	Ensure, by pulling up, that the rails are fully locked.
	There are entrapments or crushing risks when the headboard or footboard is being placed. Ensure that there are no obstacles during operation.

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* The installation and the extraction HPL headboards and footboards is done in the same way.

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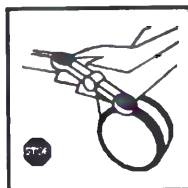
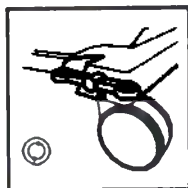
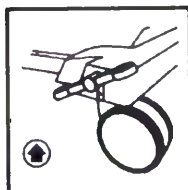
* Optionally, the bed has an unlock acoustic, which beeps when the brake pedal is not activated and the bed is still plugged.

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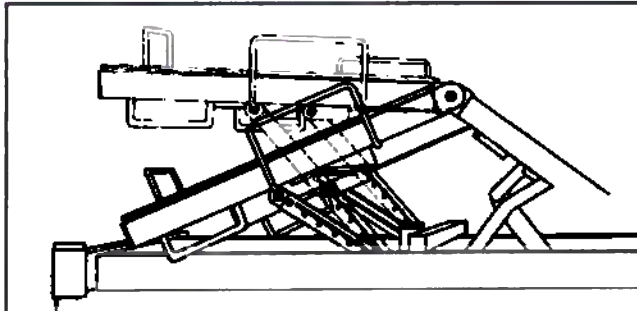
Once the legs module has been electronically risen, move the feet section to reach the desired position. The movement of the foot module can be made thanks to the rastomat fitting with 6 positions that place the module at different heights. To increase the inclination, it is sufficient to pull up the module. Make sure that both hinges are correctly positioned.

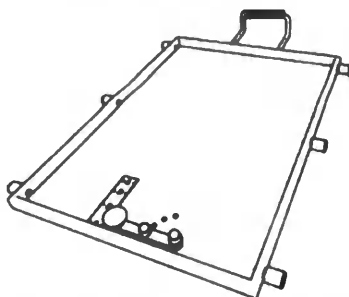
To decrease the inclination, raise the module slightly and lower it to the desired position.

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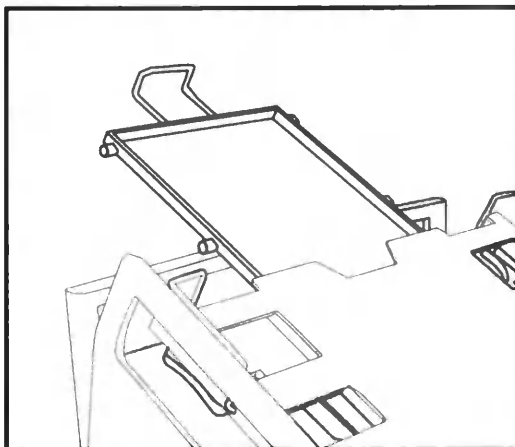
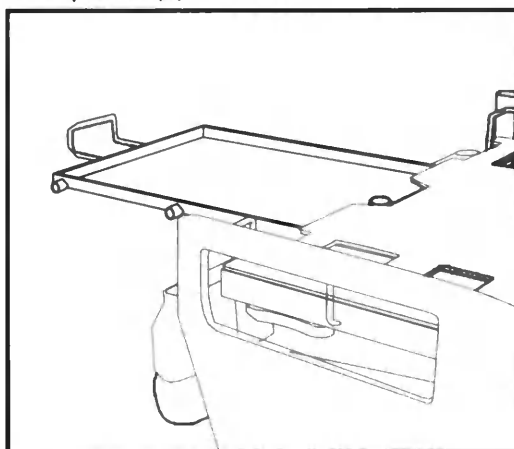
- Rise the backrest until it is over the headboard and be able to introduce it (1).

- Remove the headboard and introduce the cassette carrier without the necessity to lift the backrest avoiding to move the patient (2)..

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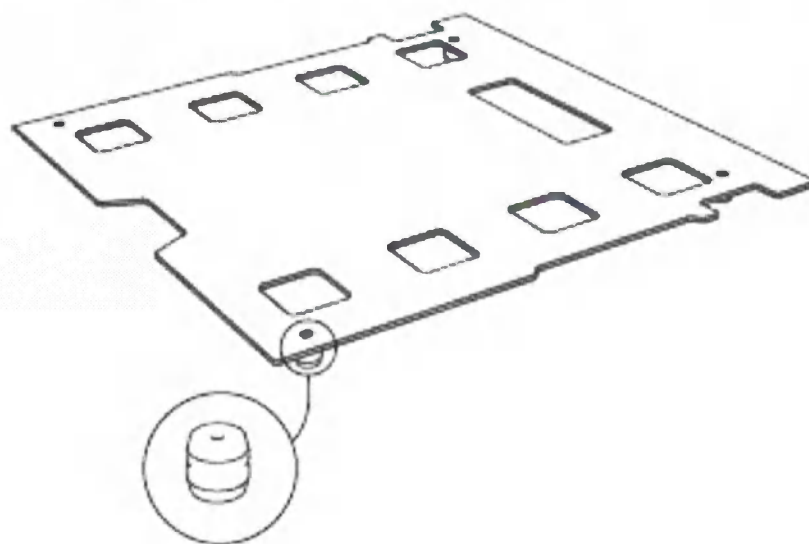
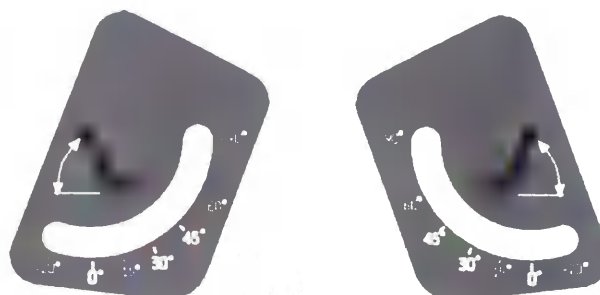
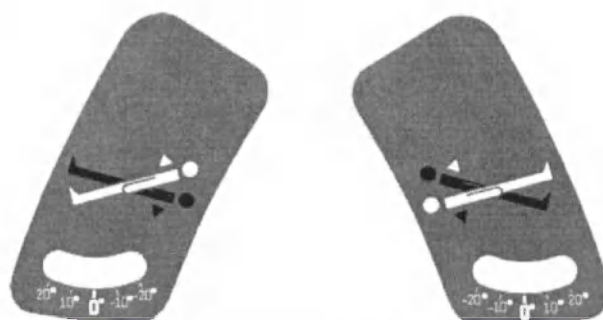
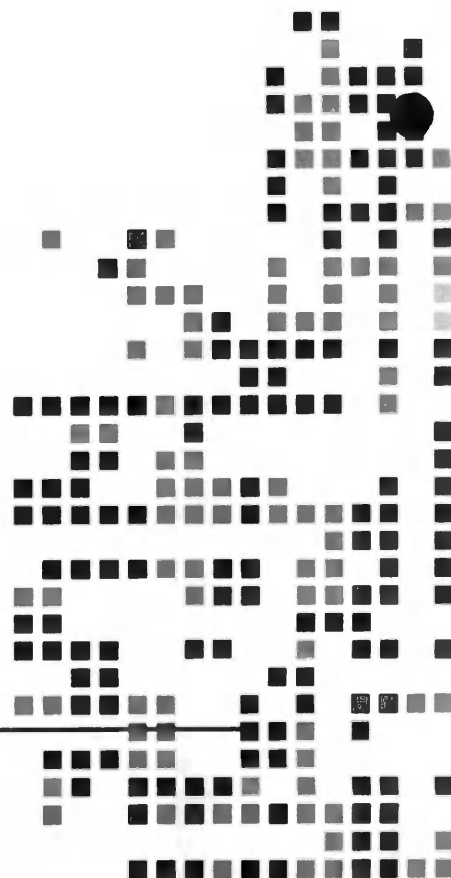
Trendelenburg and reverse trendelenburg

Located in the extremities rails, it shows the tilting degrees in the trendelenburg and reverse trendelenburg movements of the bed.

Cleaning and maintenance

Bed cleaning

*It is recommended not to use any steaming cleaning systems because an excessive humidity and temperature may damage some components.
Do not use organic dissolvent for bed cleaning.*

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Disinfection:

- Wash the surface
- Wipe the surface with a sponge or soft cloth dampened with an approved cleaning agent
- Remove residual water with a soft cloth or sponge
- Dry.

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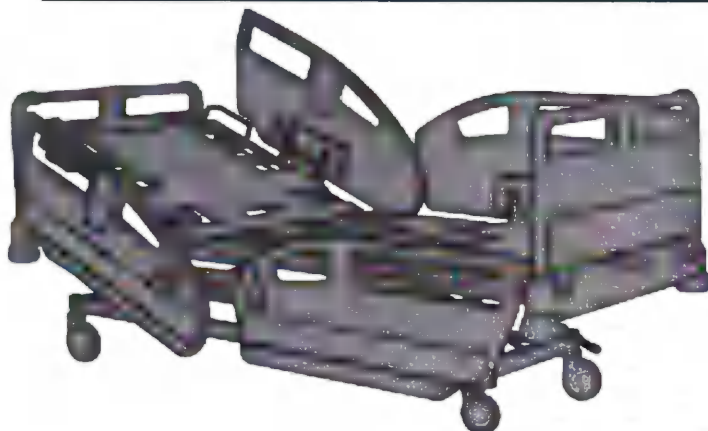
- The entire bed.

Steam cleaning

It is recommended not to use any steaming cleaning device to clean NewCare bed range. An excessive humidity may damage the bed mechanisms and components.

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To clean thick or dry stains on the floor, it may be necessary to saturate them first.

Disinfection

Dilute disinfectant and / or germicidal according to manufacturer's label instructions.

Apply the cleaner to the bed with a rag, avoiding excessive liquid, or use a surface disinfectant spray and spread the cleaner over the bed with a rag.

Allow the product to dry for the time recommended by the product manufacturer to ensure maximum effectiveness.

We recommend keeping a log of disinfection performed for each bed, documenting the date of disinfection, bed ID, time elapsed since last disinfection and disinfection methods and / or disinfectants used.

If none of the approved cleaning agents provide the required level of disinfection, use a 5.25% sodium hypochlorite solution diluted as follows: 1 part disinfectant to 100 parts water (by volume).

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All the points have to be followed making the appropriate and necessary changes during the procedures.

This guide is designed to be used along with the preventive maintenance chart. This chart is designed to keep a maintenance history and the consequent repair costs of the bed. Taking specific notes and fulfilling the maintenance can increase the bed life time.

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- Lift electrically the backrest until the maximum tilted position is reached.
- Verify that the cables that goes from the emergency RCP actioning levers to the motor are in good condition.
- Observe the correct union of the absorber to the bed and if it is in good conditions.
- Action independently one lever and check the correct lowering of the backrest as well as the good absorber functioning.
- Perform steps 1 and 4 activating in step 4 the other lever.

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- Verify the correct motor attachment to the bed.
- Lift the backrest to its highest position and afterward down it to its lowest position (flat).
- Check that no strange friction noise can be heard and no overload signal sounds whilst moving.

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Legs motor (check every two years)

- Verify the correct motor attachment to the bed.
- Lift the legs section to its highest position and afterward down it to its lowest position (flat).
- Check that no strange friction noise can be heard and no overload signal sounds whilst moving.

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Battery (check every two years)

- If necessary recharge the battery.
- Unplug the bed from mains.
- Try to perform all the beds functions and check that works correctly and there is not disfunction signs.
- Check that the battery date is not over three or four years. If it is exceeded, it is recommended to be replaced.

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Rapid problems resolution

PROBLEM	RESOLUTION
The bed beeps each time a button is pressed.	It is possible the power is not reaching the bed and it beeps to warn that the battery is discharging.
The bed has locked functions.	Some or all bed functions are locked. To unlock follow the next instructions: From the nursing control: The unlocking of the bedbase sections or bed functions is made by pressing the locking button with a lock icon (V) and following press once the button of the desired section (B, D, H, F, J). The light indicator (K) of the selected section will turn off.
Leds blinks and the bed does not work.	Wait till the leds stop blinking and check the locking situation as the preview section.
The bed is locked and nothing works.	The bed is locked and nothing works. Unplug from mains and battery. Wait 30 sec. and plug to mains again. Once the controls are recognised plug the battery.
The bed stay locked and the above procedure did not work.	In this case reset the system as follows: From the patient control: Press and hold A and C keys for a few seconds until the control beep disappears. The system resets and operates correctly again. From the nursing control: Press and hold simultaneously E and F keys for a few seconds till the system resets.
After reset the system the bed operates but when pressing any button it locks again.	If after reset the bed by pressing the controls it locks again, it is possible that there is a problem with the connexion cables and produces a short circuit in the control box and stops the bed. Replace the cable and call the technical service.
The battery does not work.	If after the battery charge, the charge indicator located in the nursing control does not shows that it is charged, check the battery condition. If it is below 17 volts, replace it for a new one.

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The batteries have a limited life that depends on the use done of them. The value from which it is impossible to recharge is 17 volts. Below this value, the battery is unrecoverable and it has to be replaced. Below 20 volts, it is already a low value and it is recommended not surpass it.

The batteries have to stay in the bed when it is plugged to mains. In case the bed is stored, it is recommended to unplug the battery and if the storing time is over 2-3 months, is better to remove the battery and leave it on a charger or use it in another bed that is in plugged, to avoid discharges.

The bed's clutch is only to lower the bed and not to lift it. The inappropriate use of the system voids the backrest motor and clutch warranty.

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The connection box has an overuse protector (approx. for 2 minutes) that is added to the reversible and not reversible thermics (80 y 95°) that has the control box as safety elements.

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When disconnecting from mains always wait a couple of minutes to connect it again and press the controls buttons to discharge the control PCB capacitors.

Product disposal

NewCare V2 and electronic accessories shouldn't be litter with other urban wastes. To protect natural resources and promote recycling, separate batteries and recycle them by local collection procedure. Other elements in NewCare V2 can be separated: all metallic elements, plastic elements, wooden elements, phenolic elements, packaging and RAEEs, to avoid environmental damages classify these elements and recycle them properly or entrust them in a residues delivery point in your city. Thus material sources recycling is promoted.

Private users can contact the establishment where they purchased NewCare V2 or competent authority to be informed about how and where they can deliver the bed to be recycled safely.

Electrical equipment must be sent to an approved electronic equipment recycling facility. Lead-acid batteries should be sent to a dedicated battery recycling company. Please comply with the EU Waste Electrical and Electronic Equipment Directive.

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The bed and / or its parts and accessories must be disposed of in accordance with the existing applicable local regulations. For the Russian Federation, it is recommended to dispose of in accordance with the rules approved by "SanPiN 2.1.7.2790-10 (Sanitary and Epidemiological Requirements for the Management of Medical Waste)" for class "A" (non-hazardous) waste, with the exception of electronic parts.

Electronic parts must be disposed of in accordance with local special regulations.

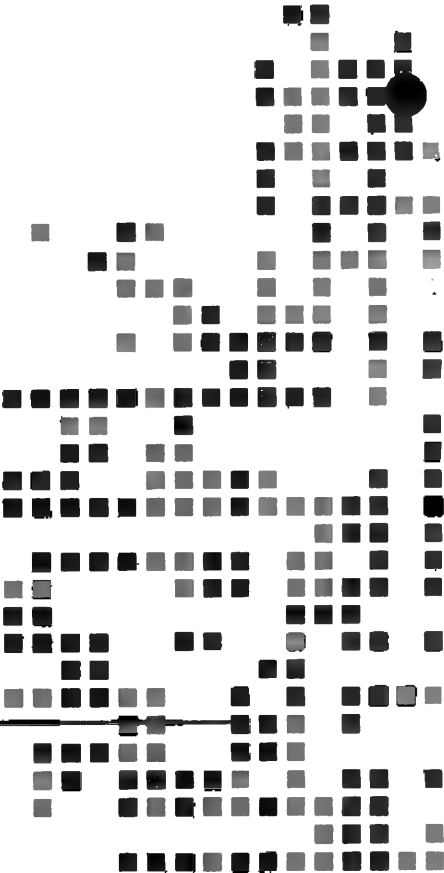
This symbol indicates separate collection of waste for electrical and electronic equipment.

Quality assurance

NewCare V2 and electronic accessories shouldn't be litter with other urban wastes. To protect natural resources and promote recycling, separate batteries and recycle them by local collection procedure. Other elements in NewCare V2 can be separated: all metallic elements, plastic elements, wooden elements, phenolic elements, packaging and RAEES, to avoid environmental damages classify these elements and recycle them properly or entrust them in a residues delivery point in your city. Thus material sources recycling is promoted. Private users can contact the establishment where they purchased NewCare V2 or competent authority to be informed about how and where they can deliver the bed to be recycled safely.

Technical specifications

Physical bed description

Measurement and weight specification	Dimension
Total length with headboard and footboard	2230 mm. $\pm$ 10 mm.
Total width with rails folded/rised	980 mm. $\pm$ 10 mm.
Minimum / Maximum bed height (without mattress) (lowest / highest position)	
• Integral wheels, diameter 150 mm:	410 mm - 815 mm $\pm$ 10 mm
• Metal wheels with individual brakes, diameter 150 mm:	395 mm - 800 mm $\pm$ 10 mm
• Metal wheels, diameter 150 mm:	415 mm - 820 mm $\pm$ 10 mm
Crane pitch height	
• Integral wheels, diameter 150 mm:	170 mm $\pm$ 10 mm
• Metal wheels with individual brakes, diameter 150 mm:	155 mm $\pm$ 10 mm
• Metal wheels, diameter 150 mm:	175 mm $\pm$ 10 mm
Bed base double regression	150 mm. $\pm$ 10 mm.
Bed section length	790 / 170 / 320 / 620 mm
Bed section width	830 mm
Backrest section maximal tilting	$65^{\circ} \pm 5^{\circ}$
Femoral section maximal tilting	$45^{\circ} \pm 2^{\circ}$
Leg section maximal tilting	$18^{\circ}/7^{\circ} \pm 2^{\circ}$
Trendelenburg / Reverse Trendelenburg	$13^{\circ} / 15^{\circ}$
Bed maximum weight	120-150 kg.
Safety weight load (SWL)	250 kg.
Maximum patient weight (MPW)	185 kg.

Rails physical description

Specifications	Dimension
Backrest rail length	1010 mm.
Legs section rail length	1040 mm.
Rails height respecting the bed base	430 mm.

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Specifications	Dimension
Mattress lengthxweight	850x2000 mm.
Mattress length	2000 mm.
Minimum mattress height	130 mm.
Maximum mattress height	210 mm.

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Condition	Range
Temperature	Between +5°C and +40°C
Humidity	
Pressure	From 500 hPa to 1060 hPa (from 0 to 2000 m. altitude)

Electric system description

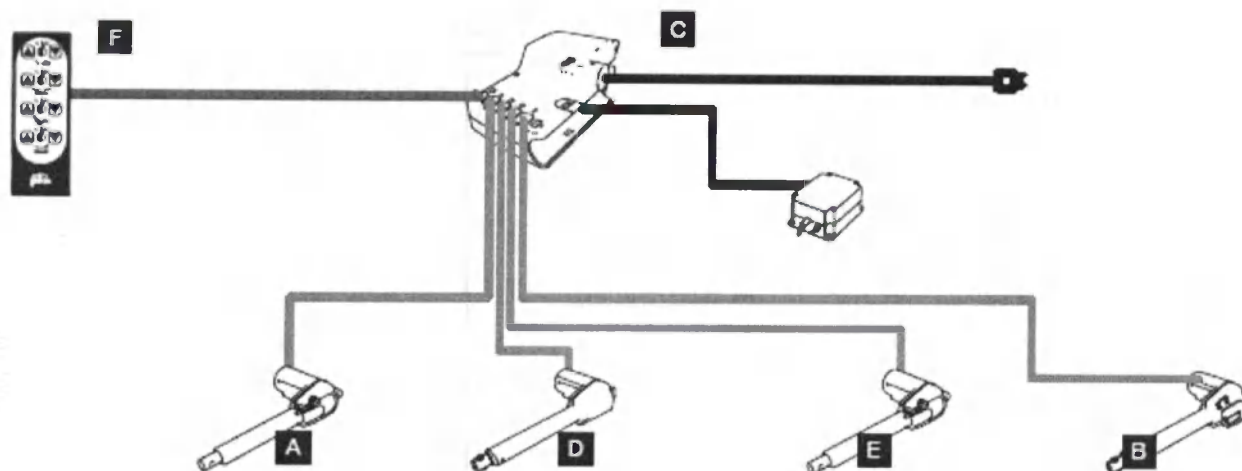
Description Specification	Specification
Electric shocks protection type (CEI 60601-1)	Class II
Electric shocks protection grade (CEI 60601-1)	Typo B
Water entry protection (CEI 60529)	IP66

Electrical diagram

NewCare V2 bed has the following configuration:

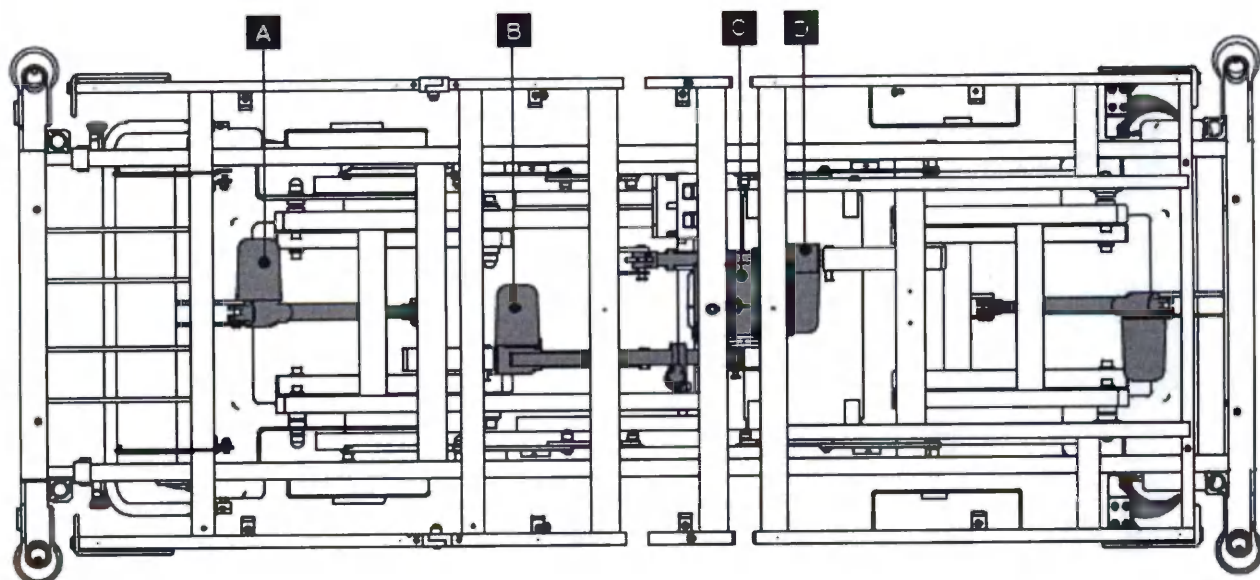
- Control box.
- Battery.
- Four lineal actuators (one for the backrest, one for the legs and two for the lifting).
- Patient control with extendable wire (depending on options).

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* This diagram is an example of a possible configuration.

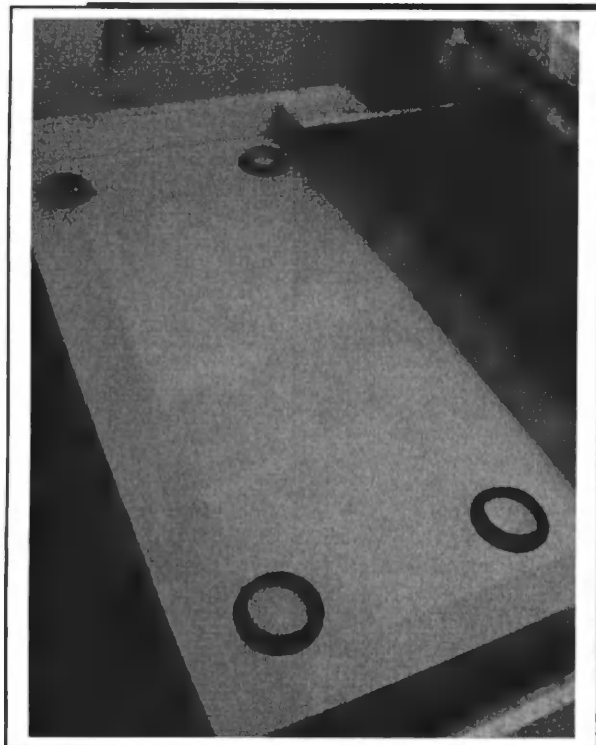
See the spare parts manual for complete diagrams of each bed configuration.

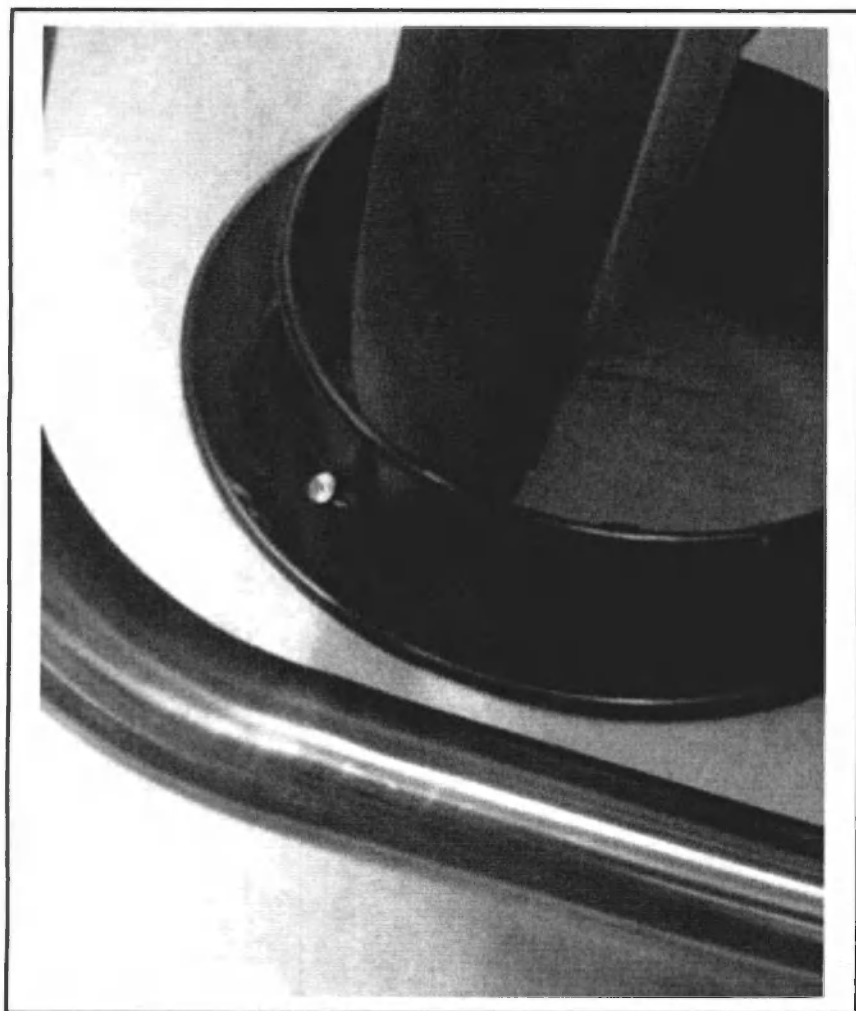
Item	Description	Item	Description
A	Bed lifting motor (femoral section)	D	Femoral section lifting motor
B	Backrest lifting motor	E	Bed lifting motor (backrest section)
C	Control box	F	Patient control

Electrical diagram

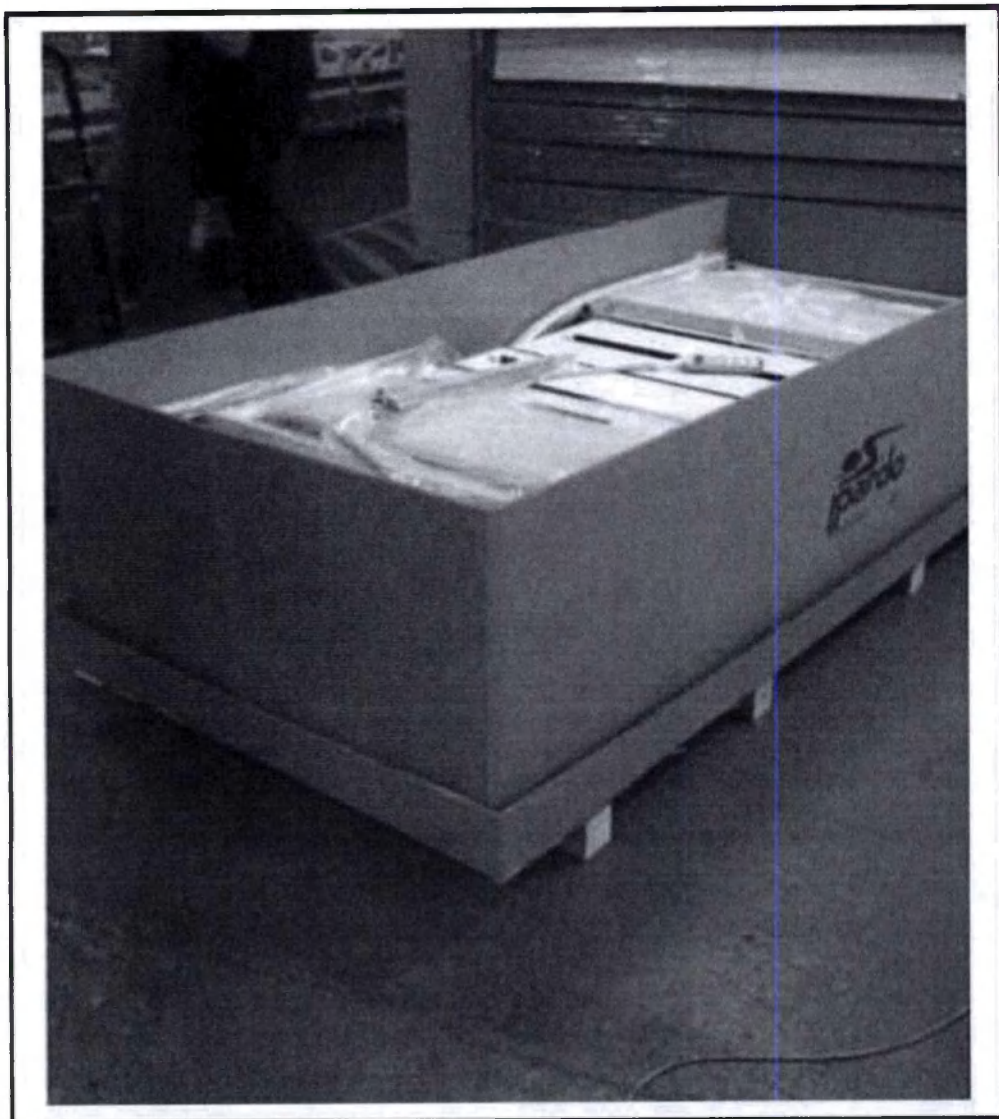
- If specified in the contract or order, the medical beds will be packed according to the model called "double lid boxes".
- The pallet is placed on the floor, one of the covers is placed on top, after closing the sashes are stapled or glued.
- Depending on the size of the bed, the mooring (fixing) wheel cylinders are placed on the first cover, as shown in the following picture:

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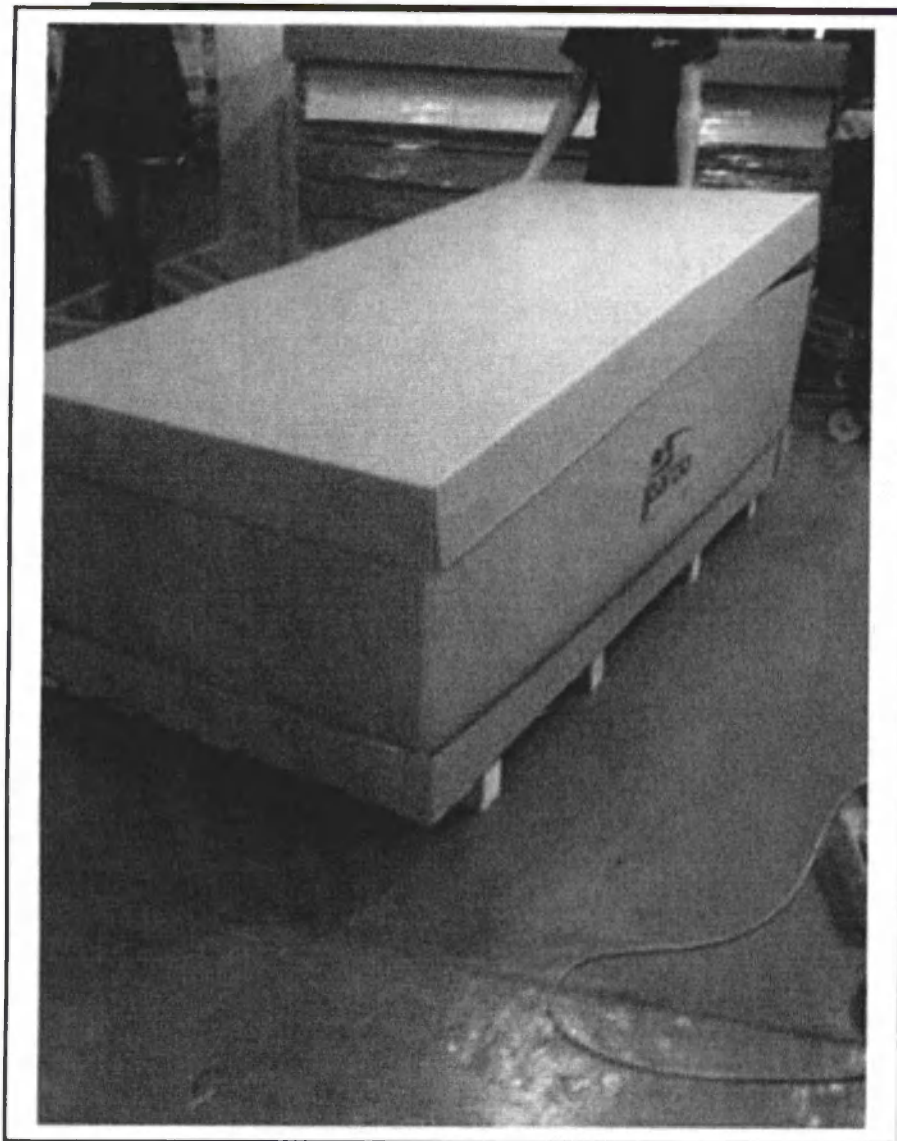
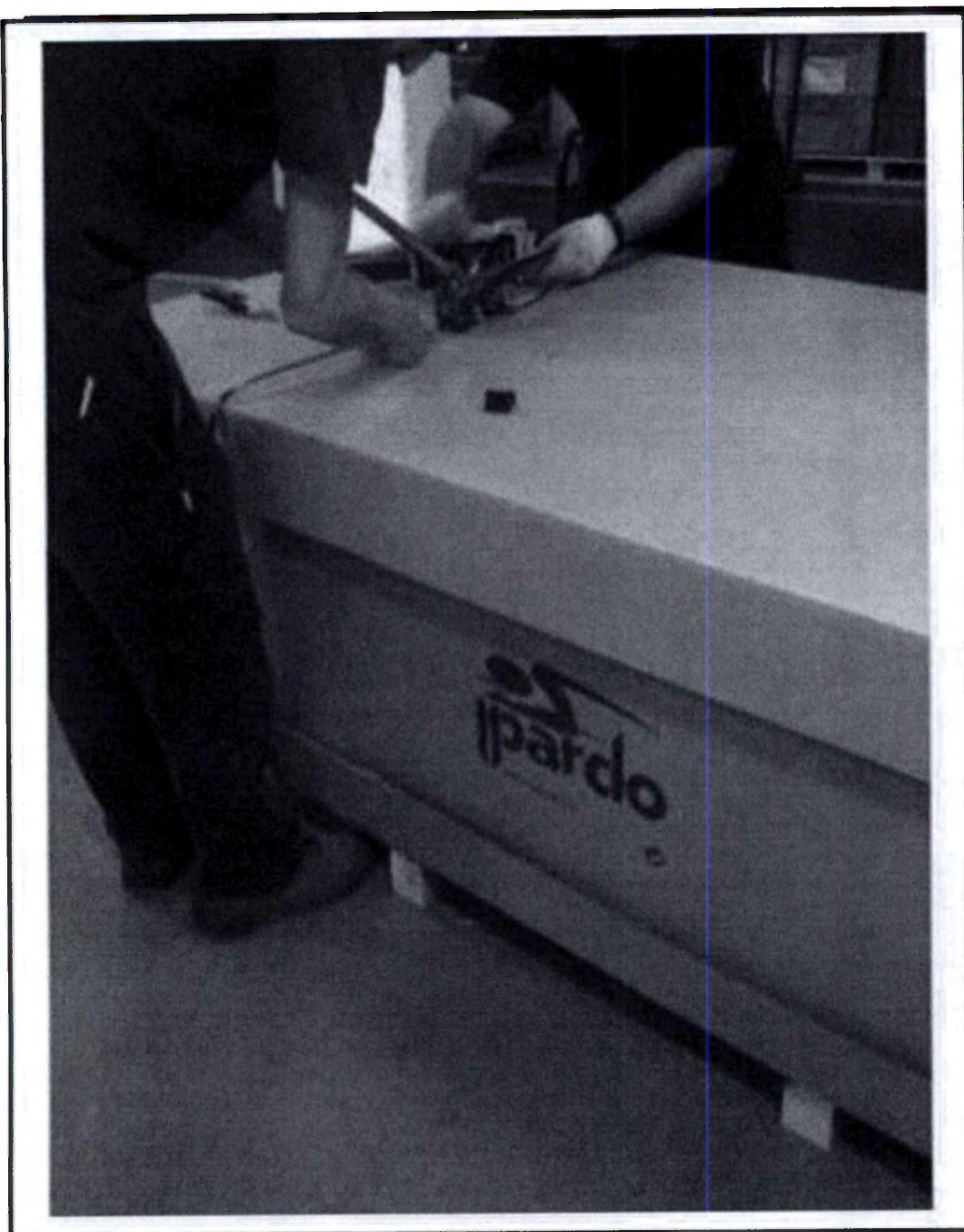
-
- The cylinders are screwed to the base.
 - The side of the box is positioned as follows:

- The cover is fitted the same size as the base.

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- For this reason, the technical service plans the specified installation as follows: clarify with the customer all data that may be of interest for the installation (for example, the specific day and time of delivery, truck access to the Center, documentation requested by the customer, etc.) ... All this information is agreed between the customer and the installer.

All additional information, where appropriate, is attached to the rest of the documentation that will be delivered to the persons assigned to this external assembly. The specified documentation will also contain at least "waybills". 3 copies of invoices are delivered: one - to the customer as acceptance of delivery, another signed by the driver as confirmation of delivery, and the third is returned by the installer as confirmation of installation.

If the exterior installation is performed by the service provider, then it must be properly assessed and approved as such. In the case of using the services of temporary employment agencies, the staff of the Pardo company

will be present (commercial, quality, SAT, etc.).

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Method of work

- Based on the data provided by the client or the conditions stipulated in the contract or order, the head of the commercial department, the person in charge of logistics, or the corresponding authorized person in accordance with the definition of the functions of the workplaces, carries out the planning of the installation, assignment of funds that he considers appropriate ...

- In addition to the relevant documentation, the personnel assigned to the installation will be issued with "delivery notes" (original + copy) of the sales item of the installation.

- Where the installation includes assembly, the necessary assembly instructions are included with the product parts that are included in the packaging step.

- The installer (installer) signs the waybill as a receipt and, in the event of an incident, sends it by mail.

Additionally and to the extent that the client accepts it, the client personally signs his consent in the same "invoice".

- A copy of the invoice remains with the client.

- The persons in charge of the installation will deliver the invoice signed by the customer.

- If an incident occurs during installation (breakdown due to transportation, accident, etc.), the persons responsible for the installation inform the technical department of the manufacturer about this, indicating, in addition to the customer's data, the nature of the malfunction, the work performed and the decision taken.

Environmental commitment

- The workers who took part in the installation must properly dispose of all waste that is generated as a result of their activities.

• Waste disposed of in urban environments (cardboard, plastic, etc.) should be sent to a container at a designated waste collection point near the installation site. If there is no such container nearby, they should send it to the company's municipal waste containers.

• If hazardous waste is generated (aerosols, contaminated rags, etc.), it should be disposed of in an appropriate location in the company's hazardous waste warehouse.

- If you have any doubts about environmental issues, you can contact an environmental technician.

"EC" declaration of conformity

59

"EC" DECLARATION OF CONFORMITY

INDUSTRIAS PARDO, S.L.
C/ Los Angeles 5, (Polígono Centrovía)
50198 La Muela, Zaragoza (Spain)
CIF: B50057322

Authorized company by the Spanish Agency of Medicines and Medical Devices (AEMPS) of the Spanish Ministry of Health nº 261 P.S.

HEREBY DECLARES THE DESIGN AND MANUFACTURING OF THE FOLLOWING PRODUCTS:

Type: MEDICAL BED
Model: NEWCARE (V2)
Serial number: See labels
Manufacturing year: See labels

This product is classified as a class I device, according to annex 93/42/CE.

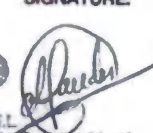
ARE IN CONFORMITY WITH THE FOLLOWING PROVISIONS:

- Directive 93/42/CE European Parliament and Council, 14th June 1993, concerning medical devices.
- Directive 2006/42/CE European Parliament and Council, 17th may 2006, concerning the machinery by amending Directive 95/16/CE.
- Directive 2006/95/CE European Parliament and Council, 12th December 2006, relative to the approximation of the legislations of the Member States on the electrical material destined to be used with certain limits of tension.
- Directive 2004/108/CE European parliament and Council, 15th December 2004, relative to the approximation of the legislations of the Member States as for electromagnetic compatibility.

STANDARDS OF REFERENCE:

- Standard UNE-EN ISO 12100 "Safety of machinery. General principles for design risk assessment and risk reduction"
- Standard UNE EN 20324, "Degrees of protection provided by enclosure, IP code". Same as IEC 60529"
- Standard UNE EN 60204 -1, "Safety of the machinery. Electrical equipment of machines. Part 1: General requirements"
- Standard UNE EN 60601-1, "Medical electrical equipment. Part 1. General requirements for safety and essential performance"
- Standard UNE EN 60601-1-2, "Medical electrical equipment. Part 2: General requirements for safety. Collateral standard: electromagnetic compatibility. Requirements and test"
- Standard EN 60601-2-52, "Medical electrical equipment. Part 2-52: Particular requirements for basic safety and essential performance of medical beds"

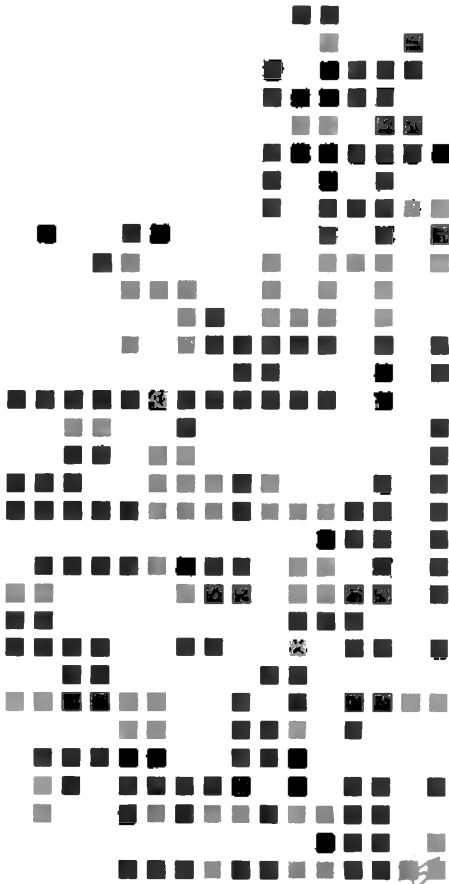
SIGNATURE:

Industrias Médicas Pardo, S.L.
NIF B-50057322 Adolfo Gaudes Martí
C/ Los Angeles 5 - 50198 LA MUELA (ZARAGOZA)
Technical Quality Manager
Tél: 076 50 20 20 Zaragoza, 04th November 2015
www.pardo.es

Preventive maintenance control list

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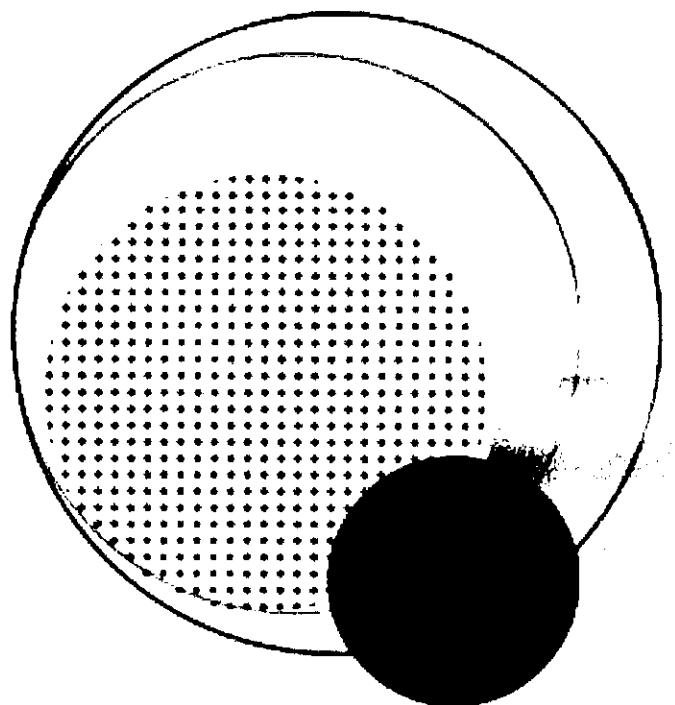
FUNCIONES								
	Date	Date	Date	Date	Date	Date	Date	Date
General Status								
Headboard and footboard								
Brake system								
Electric system								
Rails								
Wheels								
Bed base								
Sections motors								
Lifting motors								
Extensible, rastomat, RCP								
Others								

Repair cost								
-------------	--	--	--	--	--	--	--	--

Labour cost								
-------------	--	--	--	--	--	--	--	--

Total cost								
------------	--	--	--	--	--	--	--	--

Reviewed by								
-------------	--	--	--	--	--	--	--	--

Mattresses

bed models can be equipped with the following mattress models:

Model	Q-ty of layers	Foam type	Уровень защиты от пролежней	Max patient weight	Top profile	Fire resistance	Symmetry
Density 900	3	HR + Visco	I – very high	250 kg	Yes	No	Yes
Density 650 Fresh	3	HR + Visco	II – high	250 kg	No	Yes	No
Density 600	3	HR + Visco	II – high	250 kg	Yes	No	Yes
Density 475 Fresh	2	HR + Visco	II – high	230 kg	Yes	No	Yes
Density 450	2	HR + Visco	III-middle	230 kg	Yes	Yes	Yes
Density 400	2	HR + Visco	III-middle	230 kg	No	No	Yes
Density 250 Fresh	1	Biocell	IV – low	190 kg	No	No	Yes
Density 150	1	HR	IV – low	190 kg	Yes	No	Yes
Density 100	1	HR	IV – low	190 kg	Yes	No	Yes

HR - Polyurethane foam (latex free)
Visco - Memory Foam (Latex Free)
Biocell - Biocell foam BIOCEL3230 (latex free) Eco-friendly foam with high elasticity, made from natural plant products (in particular soy), renewable and environmentally friendly.

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DENSITY 100

FIRA
FEDERATION OF
INDUSTRIAL RESEARCH ASSOCIATIONS

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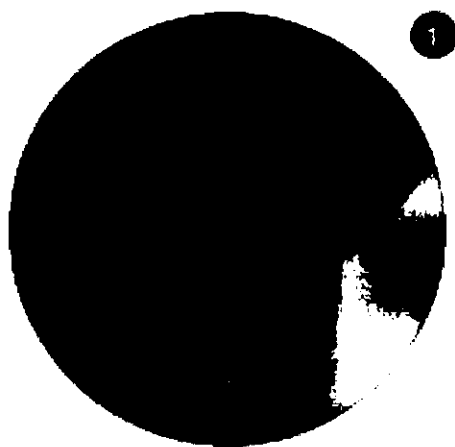
High density monolithic foam: core (25 kg/m^3), with open cell structure that improve the air exhaust through the foam, becoming the core breathable and hygienic.

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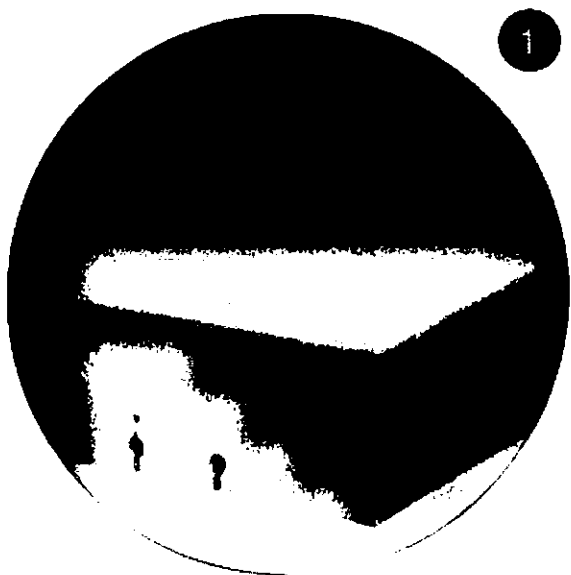
High resiliency ecologic Biocell foam (32 kg/m^3), made of vegetal origin natural products (specifically soya), renewable and non-polluting; with a height of 15 cm. thick.

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Framed lower mattress cover to maintain patient body curvature and bed articulation. High adaptability of the mattress to the bed promoting better articulation.

Totally ergonomic, keeps the body in its natural curved position, preventing problems with the back and cervical pain.

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Symmetric core allows using indistinctly the head or feet area.

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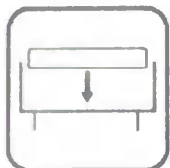
The head and foot area can be used indistinctly.

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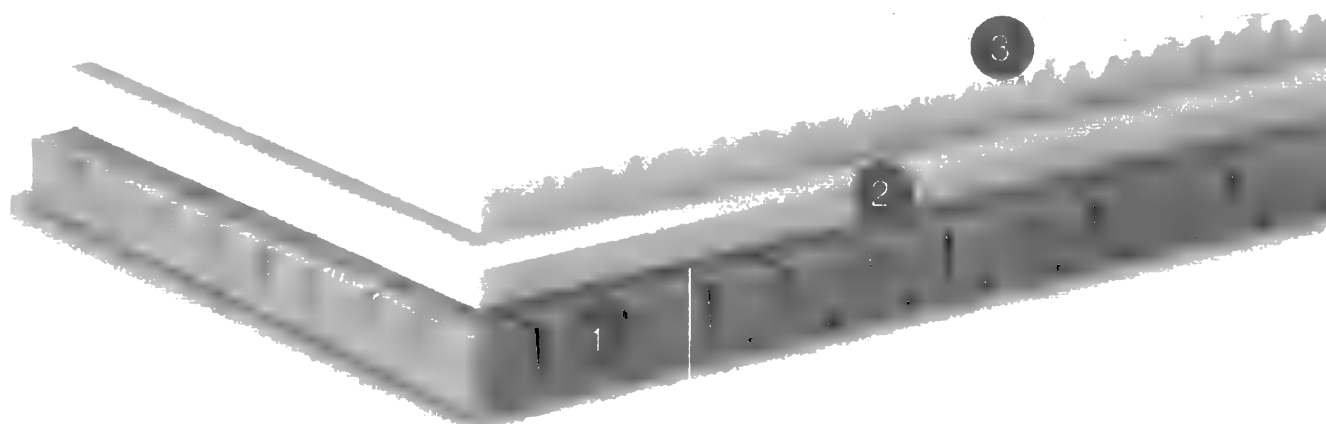
HR polyurethane foam base, 38 kg/m³, with 8 cm. of height. Its height and porosity allows a great air circulation through the foam.

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Thermosensitive viscoelastic foam 50 kg/m³, with 4 cm. height for the upper layer. When in contact with the body, it becomes softer, maintaining the same level of support, providing excellent levels of comfort.

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Version

Medical functional electric 4-section NewCare V2 Standard bed consisting of:

1. Bed frame - 1 pc.
2. Bed platform, composition:
 - removable sections of the stock made of HPL-plastic - 3 pcs.
 - fixed section of the stock made of steel - 1 pc.
3. Plastic side rails of the back section (if necessary), execution options:
 - with built-in control panel (if necessary) - 2 pcs.
 - without control panel (if necessary) - 2 pcs.
4. Plastic side rails of the leg section (if necessary) - 2 pcs.
5. Metal side rails, execution options:
 - stainless steel (if necessary) - 2 pcs.
 - painted steel (if necessary) - 2 pcs.
6. Head end, execution options:
 - removable with hand grips (if necessary) - 1 pc.
 - removable straight (if necessary) - 1 pc.
 - non-removable (if necessary) - 1 pc.
7. Foot end, execution options:
 - removable with hand grips with built-in control panel (if necessary) - 1 pc.
 - removable with hand grips without a control panel (if necessary) - 1 pc.
 - removable straight line with built-in control panel (if necessary) - 1 pc.
 - removable straight line without control panel (if necessary) - 1 pc.
8. Wheels, execution options:
 - metal 150 mm (if necessary) - 4 pcs.
 - integral 150 mm (if necessary) - 4 pcs.
 - double 150 mm (if necessary) - 4 pcs.
9. Central brake, execution options:
 - multifunctional pedal (if necessary) - 1 pc.
 - side pedals of the central brake on the wheels at the foot end (if necessary) - 2 pcs.
 - side pedals of the central brake for all wheels (if necessary) - 4 pcs.
10. Foot control pedal (if necessary) - 1 pc

11. Anti-decubitus mattress, variants:

- Density 100 (if necessary) - 1 pc.
- Density 150 (if necessary) - 1 pc.
- Density 250 fresh (if necessary) - 1 pc.
- Density 400 (if necessary) - 1 pc.
- Density 450 (if necessary) - 1 pc.
- Density 475 fresh (if necessary) - 1 pc.
- Density 600 (if necessary) - 1 pc.
- Density 650 fresh (if necessary) - 1 pc.
- Density 900 (if necessary) - 1 pc.

12. Hygienic mattress cover, variants:

- Cugar (if necessary) - 1 pc.
- Select (if necessary) - 1 pc.

13. Lever for emergency lowering of the back section - 2 pcs.

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