I ROGER HARRY NARRINGTON NOTARY PUBLIC
CERTIFY THAT I TODAY SAW KAMALJIT HUNJAN
EXECUTE THIS DOCUMENT HE HAVING PROVIDED
EVIDENCE OF AUTHORITY AND IDENTITY TO ME

R. Narrington
Notary Public 5/5/2022
Roger Narrington
Notary Public
12 Broadway
Chinham
Essex RM13 8YL
Tel: 01708 552804
Email: r.narrington@sackvilles.co.uk

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KeyMed (Medical & Industrial Equipment) Ltd

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KeyMed House, Stock Road, Southend-on-sea, Essex, SS2 5QH

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Instructions for Use:


“WM-NP3 Mobile Workstation and Accessories”

OLYMPUS

INSTRUCTIONS FOR USE

MOBILE WORKSTATION

WM-NP3

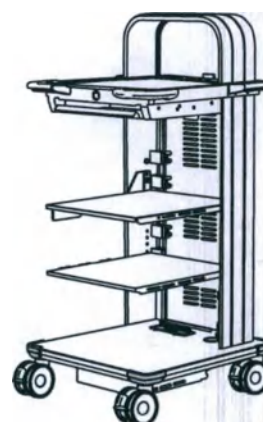K10029987
K10035360
K10035361
K10035362
K10035363

K10035364
K10035365
K10035366
K10035367
K10035368

K10035370
K10035371
K10035373
K10037157
K10037158

K10037159

Caution: Federal (USA) law restricts this device to sale by or on the order of a physician.

Contents

1	Using this document	5
1.1	User instructions	5
1.2	Conventions throughout this document	5
1.3	Safety-related information in this document.....	6
2	Safety information	7
2.1	Intended purpose	7
2.2	Contact	8
2.3	General safety notes.....	9
3	Product description	11
3.1	Standard sets.....	11
3.2	Optional accessories	12
3.3	Scope of delivery	14
3.4	Labels and symbols.....	15
3.5	Power cable	20
3.6	Warranty	21
4	Before use	22
4.1	Handling pre-installed features	22
4.2	Removing and attaching the rear panel	30
4.3	Installing MAJ-2159 Scope Pole Kit	32
4.4	Installing keyboards	43
4.5	Installing medical electrical equipment	44
5	Maneuvering.....	51
5.1	Specific safety notes for maneuvering	51
5.2	Preparing for maneuvering	51
5.3	Maneuvering the WM-NP3 Mobile Workstation	52
6	Use	54
6.1	Specific safety notes for use	54
6.2	Using the WM-NP3 Mobile Workstation.....	56

7	Reprocessing	58
7.1	Safety notes for reprocessing	58
7.2	Using reprocessing agents	58
7.3	Cleaning and disinfection.....	60
8	Maintenance	63
9	Storage and disposal.....	65
9.1	Storage	65
9.2	Disposal	65
10	Technical data	66
10.1	Power cable	66
10.2	Socket adaptor cable	67
10.3	Protective earth terminal	67
10.4	Winding.....	68
10.5	Primary winding	68
10.6	Secondary winding	68
10.7	Potential equalization terminals	69
10.8	Block diagram	70
10.9	Safety standards for medical electric systems.....	71
10.10	EMC Information.....	73
10.11	Specifications.....	76
11	Customer spare parts.....	80
11.1	WM-NP3 Mobile Workstation spares.....	80
11.2	Scope pole kit spares	92
11.3	Sliding keyboard tray spares.....	95
11.4	Rear panel spares	96
12	Troubleshooting.....	99
13	End of life.....	100
14	Glossary.....	101
15	Index.....	103

1 Using this document

This chapter provides general information on how to use these instructions for use.

1.1 User instructions

Before use, thoroughly read these instructions for use and the instructions for use of all other items that will be used during the procedure.

If the required instructions for use are missing, immediately contact an Olympus representative.

Keep the instructions for use in a safe, accessible location.

1.2 Conventions throughout this document

	This is the safety alert symbol. It is used to alert the user to potential physical injury hazards. Observe all safety messages that follow this symbol to avoid possible injury.
	This symbol indicates additional helpful information.
1. 2. 3.	A numeration indicates a sequence of actions.
•	Bullet points indicate individual actions that do not have to be performed in a specific order.
►	Arrows indicate the results of actions.
- - -	Dashes indicate the listing of data, options or objects.
1) 2) 3)	Numbers with right parenthesis name elements in illustrations

1.3 Safety-related information in this document

General safety notes

General safety notes address risks that apply to the general use and handling of the product.

Specific safety notes

Specific safety notes address safety issues that are specific to the information and actions described in a chapter or section. Specific safety notes are located at the beginning of the corresponding chapter or section.

Safety messages

Safety messages focus on unexpected incidents that can occur during the intended purpose. Safety messages are placed in the context of such risks and actions. Safety messages are highlighted by the following signal words and symbols.

WARNING

Indicates a potentially hazardous situation which, if not avoided, can result in death or serious injury.

CAUTION

Indicates a potentially hazardous situation which, if not avoided, can result in minor or moderate injury.

NOTICE

Indicates a property damage message.

2 Safety information

This chapter provides general safety information for using the product. The safety information includes the following.

- Intended purpose
- Manufacturer information and contact details
- General safety notes

2.1 Intended purpose

2.1.1 Intended use and indications

The WM-NP3 Mobile Workstation Narrow Platform 3 is intended to assist endoscopic procedures using Olympus medical devices and their accessories within medical facilities by providing a stable mains power supply and a platform for the efficient use and handling of the devices during endoscopic diagnostic and therapeutic procedures. The WM-NP3 workstation can also allow visualisation during the procedure when connected to a monitor. Do not use the device for any other purpose.

2.1.2 Contraindications

There are no known contraindications.

2.1.3 Patient target group

Not dedicated to any patient population.

2.1.4 Intended user

Medical use

Only for use by a qualified physician or healthcare professional in a professional healthcare facility.

Reprocessing

Reprocessing of reusable products may only be performed by qualified hygiene personnel.

Repair

Repair of the product may only be performed by trained qualified servicing personnel that has been authorized by Olympus. Otherwise, Olympus cannot be held responsible for the safety and performance of the product.

2.1.5 Clinical benefit

The clinical benefit of the device lies in its functionality as described by its intended purpose.

2.2 Contact

2.2.1 Manufacturer

KeyMed (Medical and Industrial Equipment) LTD
KeyMed House, Stock Road,
Southend-on-Sea, Essex SS2 5QH
United Kingdom
www.olympus.com

2.2.2 Incident reporting

European Union

If a serious incident occurs with the device, report it to the manufacturer and the relevant national authority.

Eurasian Economic Union

If an undesirable event occurs with the device, report it to the manufacturer or the authorized representative.

Other regions

If a serious incident occurs with the device, report it to the manufacturer and/or the relevant national authority according to national legislation.

2.3 General safety notes

The following safety information addresses points that apply to the general handling and use of the product. This information is to be supplemented by the information given in each chapter of this document.

CAUTION

Risk of Injury

Exceeding the maximum load capacity can cause structural failure and instability, resulting in injury. Refer to the maximum load capacity as specified in “10.11 Specifications” on page 76.

MAJ-1661 LCD Arm - Adjustable

The WM-NP3 Mobile Workstation URO set has additional weight within its base for the use of the MAJ-1661 LCD Arm - Adjustable, as this accessory requires additional ballast weight. Refer to the MAJ-1661 instructions for use for further information on product assembly and operation. The MAJ-1661 LCD Arm - Adjustable is only compatible with the URO set. Do not use the MAJ-1661 LCD Arm - Adjustable with the Basic or GI Standard sets, as these are not fitted to carry the additional weight.

NOTICE

Risk of equipment damage

Incorrect work practices when handling this product can result in equipment damage. Observe the following.

- Only use the product in the environment specified in “Specifications” on page 76.
- Only use products specified in “Standard sets” on page 11 and “Optional accessories” on page 12 with this product.
- Do not exceed the maximum load capacity as specified in “10.11 Specifications” on page 76.

Reprocessing of reusable products

Reusable product must be reprocessed according to the applicable national and local standards and regulations.

3 Product description

This chapter provides information about the product and associated accessories.

3.1 Standard sets

Contact Olympus if any parts are damaged or missing. The WM-NP3 Mobile Workstation is available in 3 standard sets. Refer to the following table to check that all items in the set ordered are present.

Product	Standard sets		
	GI Standard	Basic Set	URO Set
	UK, Europe, USA, Japan, Rest of the World	UK, Europe, USA, Japan, Rest of the World	Europe, USA, Rest of the World
WM-NP3 Mobile Workstation	1	1	1
Workstation fitting kit	1	1	1
Separation transformer (market specific)	1	1	1
MAJ-2160 Irrigation Tube Holder	1	1	1
MAJ-2146 Sliding Keyboard Tray	1	-	1
MAJ-2159 Scope Pole Kit	1	-	-
Endoscope pole fitting kit	1	-	-
MAJ-2216 LCD Monitor Arm	1	-	-
MAJ-2211 Side Handles	1	-	1
IEC lead set	1	1	1
Socket adaptor cable (USA)	1	1	1
Socket adaptor cable (Japan)	1	1	-
Instructions for use	1	1	1

3.2 Optional accessories

The following products can be used as optional accessories.

NOTICE

Risk of equipment damage


The fitment of some accessories may not be possible when other accessories are fitted and can result in equipment damage. For any questions, please contact the Endoscope Service Center, a service center designated by Olympus or an Olympus branch or sales office.

Code	Description
K10036766	Socket adaptor cable (USA)
K10036760	Socket adaptor cable (Japan)
MAJ-1635	EUS connector holder
MAJ-1639	Single CO ₂ holder (up to Ø140 mm)
MAJ-1642	IV Pole
MAJ-1650	Single CO ₂ holder (Ø140 mm to Ø205 mm)
MAJ-1653	Double CO ₂ holder
MAJ-1654	Equipotential terminal strip
MAJ-1657	Keyboard arm - side installed
MAJ-1661	LCD Adjustable Arm (only compatible with the URO set)
MAJ-1663	EUS Arm mount kit
MAJ-1665	Camera head holder
MAJ-1685	IV Bag Hanger
MAJ-1938	Holder for HC-RP / UPD-3
MAJ-1985	Cylinder Hose with Switch Over Valve (PIN) / UHI-4
MAJ-1986	Cylinder Hose with Switch Over Valve (DIN) / UHI-4
MAJ-2145	ENT Scope Hanger
MAJ-2146	Sliding Keyboard Tray
MAJ-2149	Dual Monitor Arm
MAJ-2154	Consumable Stowage Holder
MAJ-2158	Suction Jar Holder
MAJ-2159	Scope Pole Kit
MAJ-2160	Irrigation Tube Holder

Code	Description
MAJ-2165	Drawer Unit
MAJ-2166	Sterile Water Holder
MAJ-2167	Nurses Control Arm
MAJ-2173	Side Shelf
MAJ-2211	Side Handles
MAJ-2216	LCD Monitor Arm 6.5-12 kg
MAJ-2217	LCD Monitor Arm 12-14 kg

3.3 Scope of delivery

The following section describes the scope of delivery of the GI standard set. Refer to the "Standard sets" on page 11 for details on other sets.

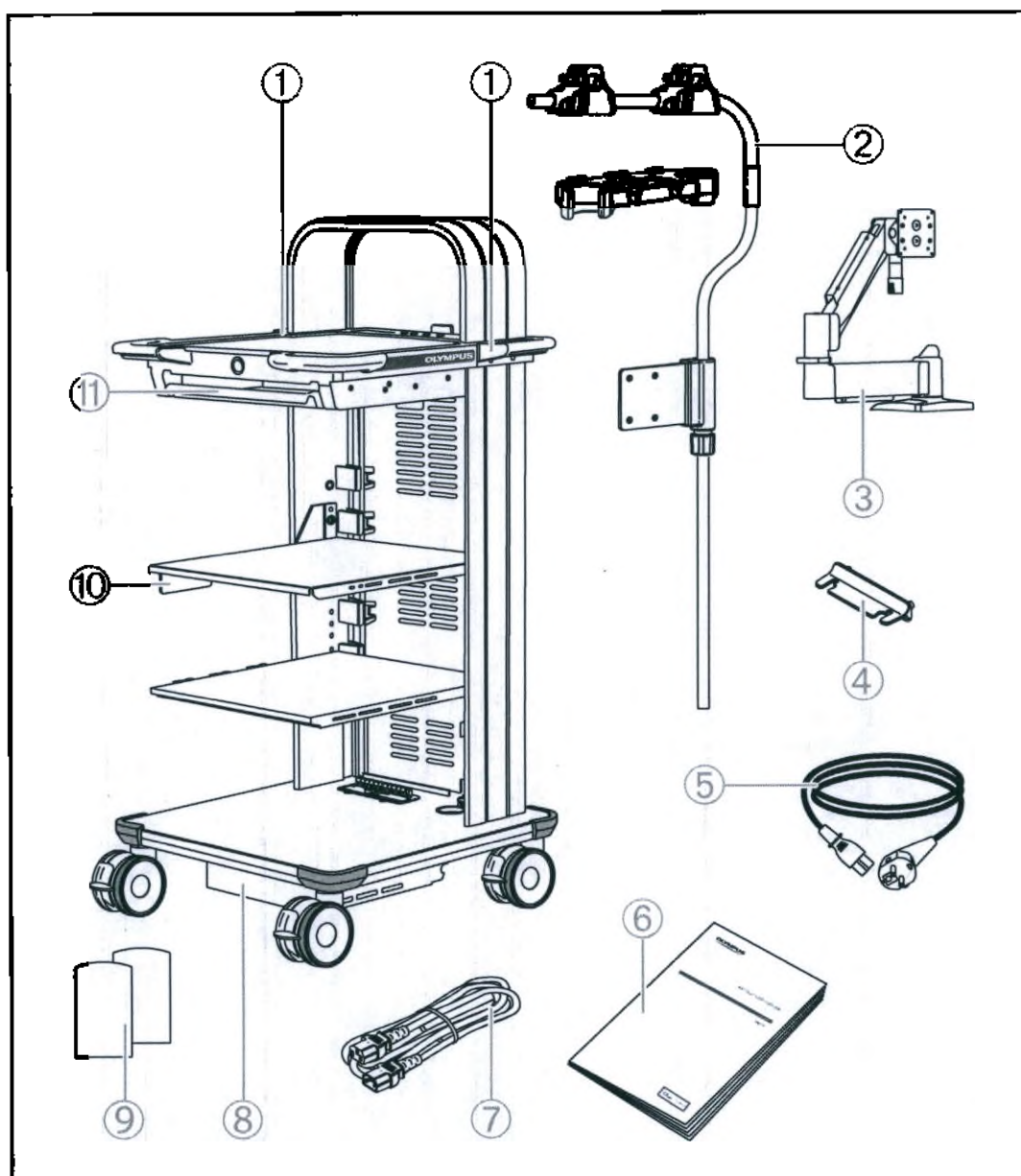- 1) MAJ-2211 Side Handles (fitted)
- 2) MAJ-2159 Scope Pole Kit
- 3) MAJ-2216 LCD Monitor Arm
- 4) MAJ-2160 Irrigation Tube Holder (fitted)
- 5) Socket adaptor cable (Japan)


- Socket adaptor cable (USA)
- 6) WM-NP3 Mobile Workstation instructions for use
- 7) IEC lead set, containing:
 - 3 x cable universal IEC 320 0.8 m
 - 5 x cable universal IEC 320 1.5 m
 - 1 x cable universal IEC 320 M-FEM 2.5 m (in separate bag)
- 8) Separation Transformer (fitted)
- 9) Fitting kit for WM-NP3 Mobile Workstation, containing:
 - 1 x 4 mm hex key
 - 1 x 6 mm hex key
 - 12 x cable ties
- Fitting kit for MAJ-2159 Scope Pole Kit, containing:
 - 4 x M6 T-nuts
 - 2 x M6 button-head screws
 - 4 x M6 countersunk screws
 - 2 x M6 crinkle washer
 - 2 x M6 plain washer
 - 1 x 2.5 mm hex key
 - 1 x 4 mm hex key
 - 1 x 5 mm hex key
- 10) Water bottle bracket (fitted)
- 11) MAJ-2146 Sliding Keyboard Tray (fitted)

3.4 Labels and symbols

This section provides information on the labels and symbol on the product and on the packaging of the product.


3.4.1 Unique Device Identification

Unique Device Identification is a label required by some countries' regulations regarding identification of medical device.

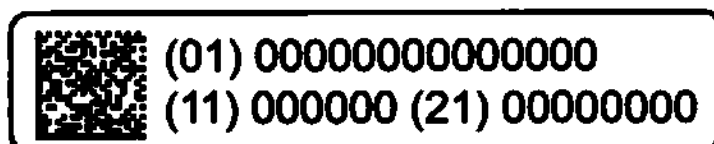The following information is coded in the 2-dimensional barcode (GS1 Data Matrix).

- (01) 14-digit GS1 Global Trade Item number.
- (11) 6-digit date of manufacture.
- (21) 8-digit serial number.

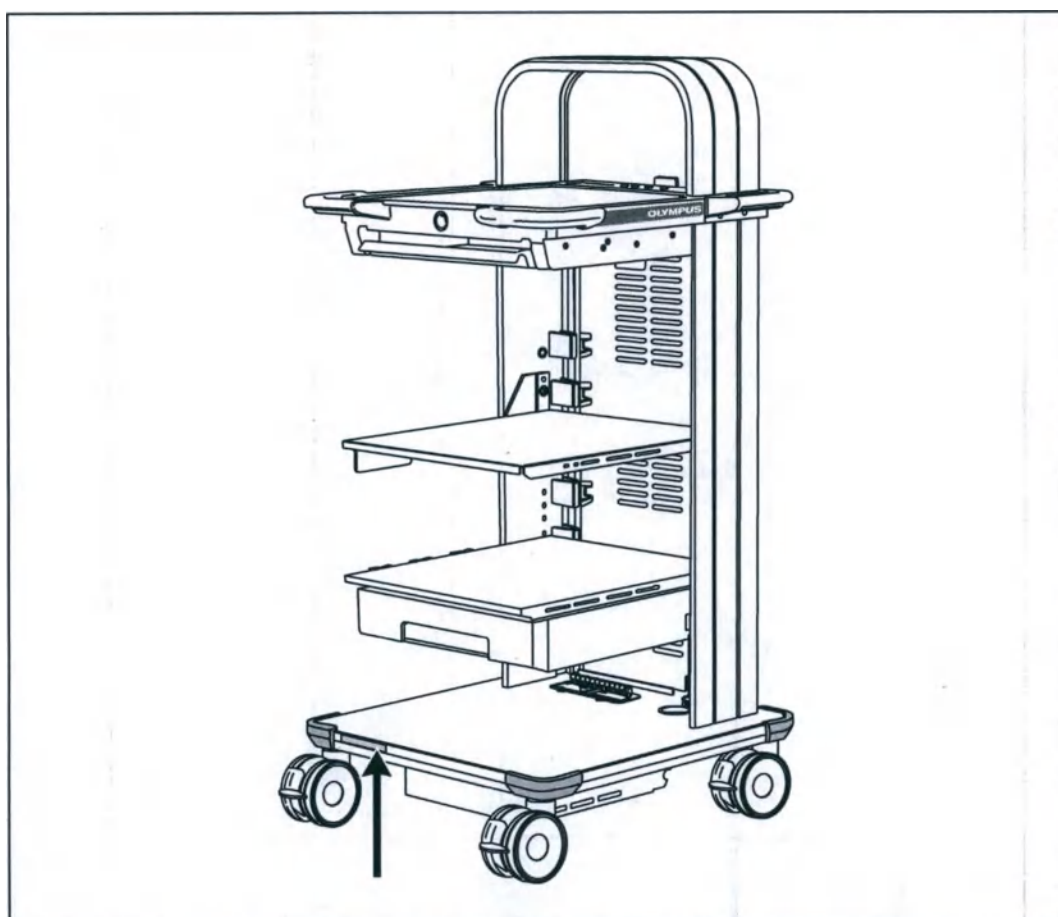
The UDI can be found on the base panel at the front of the WM-NP3 Mobile Workstation.

3.4.2 Symbols on the workstation


Symbol	Description	Source
	Refer to instructions for use	ISO 7010-M002
	Mains power ON / OFF switch	IEC 60417-5010
 Medical Device	MDR EU regulatory requirement (Class I)	Manufacturer
	Do not push the workstation when the brakes are applied	ISO 7010-P017
 311 kg 686 lbs	Max laden weight	Manufacturer
	Unit complies with China RoHS 2 requirements	China RoHS 2 – Management methods for the Restriction of the Use of Hazardous Substances in Electrical and Electronic Products, January 2016
	Stability warning	Manufacturer
	Unit complies with METI-PSE regulations – only applicable for Japan	METI-PSE regulations (Electrical Appliances and Materials Safety Act) Japan
	Compliance to the Korean Radio Waves Act (Clause 3, Article 58-2)	Korean Radio Waves Act (Clause 3, Article 58-2)
	CE marking indicating that the device is in conformity with the applicable requirements set out in applicable European Union harmonization legislation.	MDR (EU) 2017/745

3.4.3 Packaging symbols

Symbol	Description	Source
	Consult instructions for use	ISO 15223-1
	eIfU - Instructions for use available online	ISO 15223-1
	Catalogue number	ISO 15223-1
	Keep dry	ISO 15223-1
	Manufacturer	ISO 15223-1
	Keep away from sunlight	ISO 15223-1
	Stacking limit by number	ISO 7000
	Fragile, handle with care	ISO 15223-1
	This way up	ISO 7000
	Lot number	ISO 15223-1
	Temperature limitations	ISO 15223-1
	Serial number	ISO 15223-1
	Humidity limitation	ISO 15223-1
	Atmospheric pressure limitation	ISO 15223-1
	MDR EU regulatory requirement (Class I)	Manufacturer

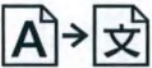	Unit complies with China RoHS 2 requirements	China RoHS 2 – Management methods for the Restriction of the Use of Hazardous Substances in Electrical and Electronic Products, January 2016
	Translation	ISO 7000
	Authorized representative in the European Community / European Union	ISO 15223-1
	UKCA marking indicating that the device is in conformity with the applicable requirements set out in the applicable UK legislation.	MHRA
	Importer	ISO 7000

3.4.4 Symbols on the separation transformer

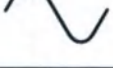Symbol	Description	Source
	Refer to instructions for use	ISO 7010-M002
	Warning to observe rules regarding the use of multiple socket outlets	ISO 7010-W001
	Equipotentiality terminal	IEC 60417-5021
	Input	IEC 60417-5034
	Output	IEC 60417-5035
	Alternating current	IEC 60417-5032
	Weight	ISO 7000-430
	Remote switch	Manufacturer

Product description

Symbol	Description	Source
	Non-inherently short circuit proof separation transformer	IEC 60417-5220
	According to European Directive 2012/19/EU on Waste Electrical and Electronic Equipment, this symbol indicates that the product must not be disposed off as unsorted municipal waste, but should be collected separately.	EN 50419
	Indicates the unit has been independently assessed for compliance with UL and CSA standards as a Recognized Component of The WM-NP3 Mobile Workstation (MAJ-2156 only).	EN 60601
	Unit complies with METI-PSE regulations – only applicable for Japan	METI-PSE regulations (Electrical Appliances and Materials Safety Act) Japan

3.5 Power cable

The WM-NP3 Mobile Workstation includes a power cable appropriate to the model of separation transformer fitted and the intended market area, as listed in "Technical data" on page 66.

- MAJ-2155 (Japan): Always use the supplied power cable, IEC, and Socket Adapter cable. Use of any other cables may cause a malfunction of the instrument or burns on the cables. Additionally, the supplied cables are for exclusive use with this instrument and must not be used with other devices.
- MAJ-2156/ MAJ-2156-15A (USA): The power cable on this unit is fitted with a 'hospital grade' plug. Grounding reliability is achieved when this equipment is connected to an equivalent receptacle marked 'hospital grade'.
- MAJ-2157 (Rest of the World): The power cable on this unit is terminated in flying leads, for a suitable plug to be fitted according to the area of intended use. Only suitably qualified personnel should fit the power plug to this cable. The cores of the power cable are colored according to the following code.
 - Brown: live ('hot')
 - Blue: neutral ('cold')
 - Green/yellow: earth ('ground')

WARNING

Risk of electric shock

Connecting the power cable to a mains electricity without protective earth can result in electric shock.

Always connect the power cable to a mains electricity with protective earth.

3.6 Warranty

Any warranty claims towards Olympus are forfeited if the user or unauthorized persons attempt repair or modification of the product. No warranty is provided for any damage due to misuse of the product.

4 Before use

This chapter provides information for installation and connection of the WM-NP3 Mobile Workstation and its accessories.


4.1 Handling pre-installed features

4.1.1 Braked castors

The WM-NP3 Mobile Workstation has braked castors. To maintain the antistatic effectiveness of the braked castors, observe the following.

- Avoid contact of the braked castors with oil.
- Avoid a build up of floor cleaning wax and polish.

Using the castors

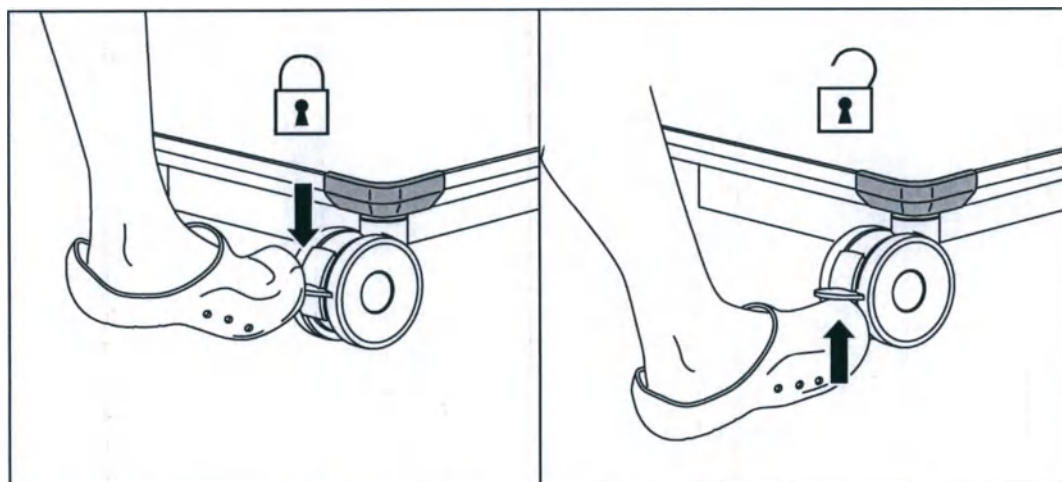
- Lock and release the braked castors as shown.
- Ensure all braked castors are released before moving the WM-NP3 Mobile Workstation.
- Always lock braked castors when the WM-NP3 Mobile Workstation is not being maneuvered.


4.1.2 Side handles

When fitted (GI and URO set), the side handles can be moved up or down to the preferred position. The side handles can also be removed if required.

Adjusting the side handles

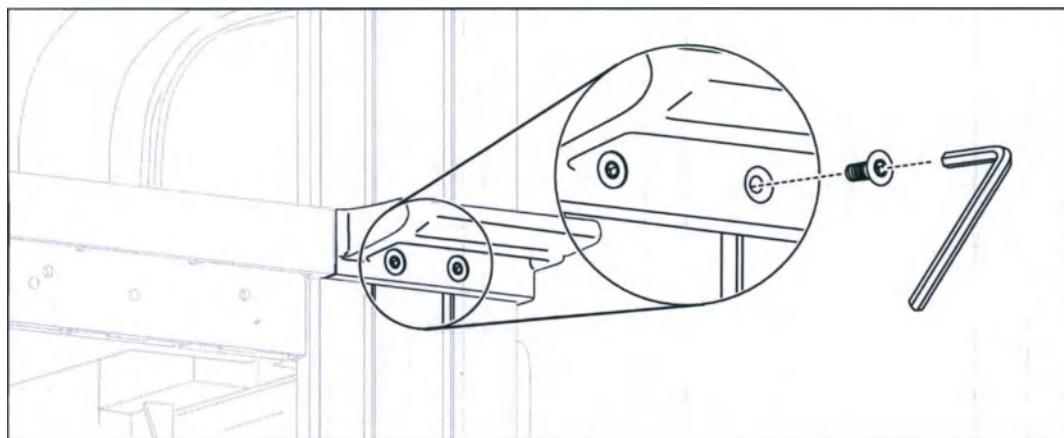1. Loosen the 2 M6 countersunk screws with the 4 mm hex key.

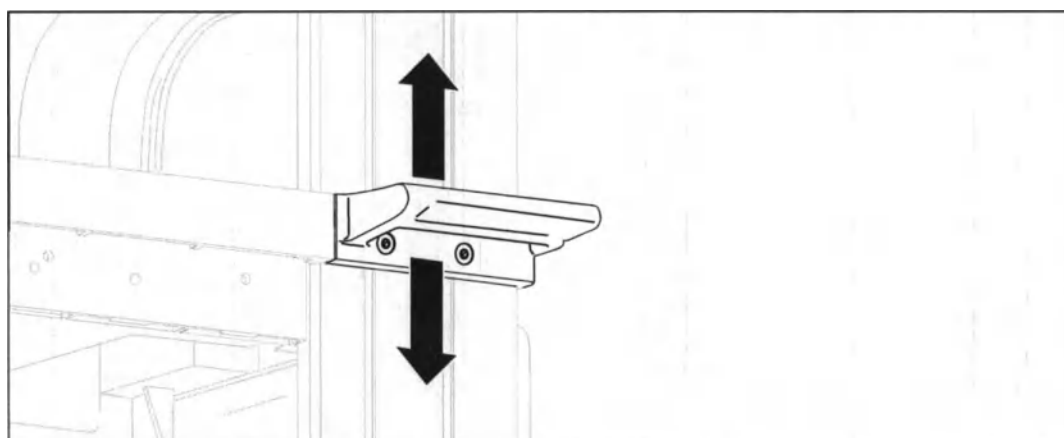2. Move the side handle to the required height.
3. Tighten the M6 countersunk screws with the 4 mm hex key.


4.1.3 Shelves

The WM-NP3 Mobile Workstation is delivered with the following shelves.

- Top tray
- 2 intermediate shelves
- Base panel


Maximum size of equipment for shelves

Refer to the maximum size of equipment (including height) when placing medical electrical equipment on the top tray, the intermediate shelves and the base panel.

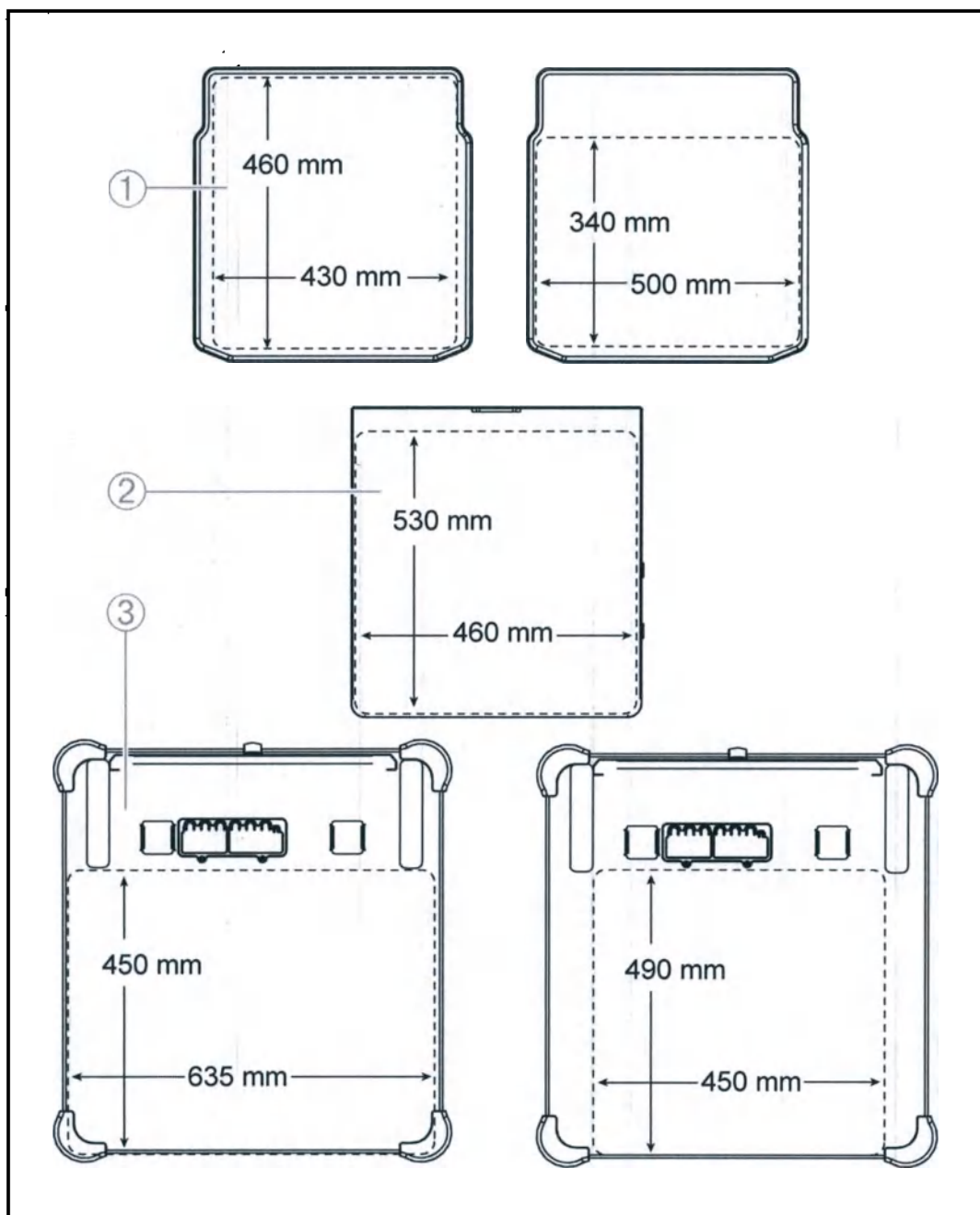- 1) Maximum size of equipment for the top tray
- 2) Maximum size of equipment for the 2 intermediate shelves
- 3) Maximum size of equipment for the base panel

Shelf positions on delivery

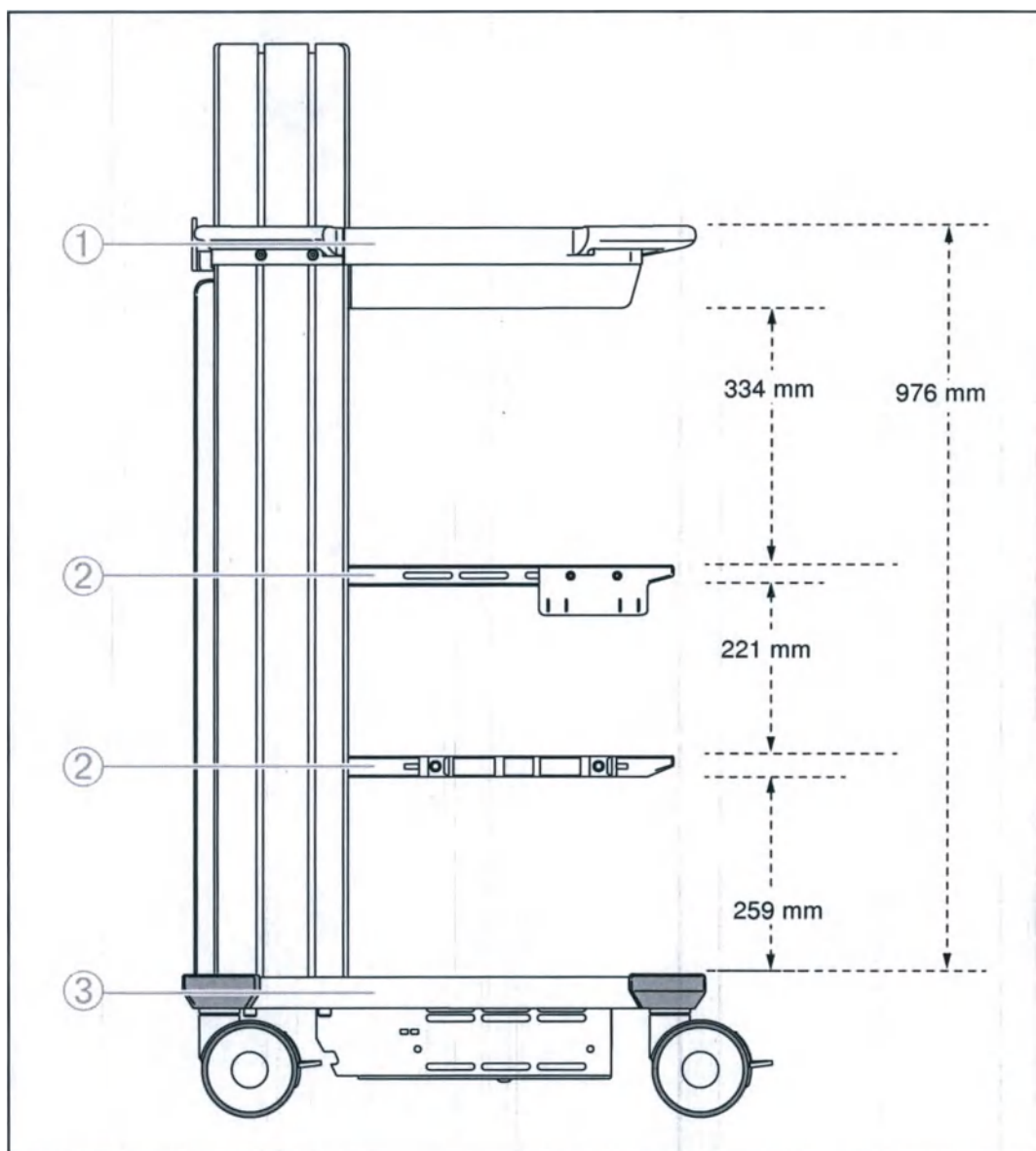
- 1) Top shelf
- 2) Intermediate shelves
- 3) Base panel

The WM-NP3 Mobile Workstation is delivered with shelves fitted in the positions shown above.

Repositioning the intermediate shelves

The columns contain holes for repositioning the intermediate shelves. The intermediate shelves are fitted with brackets that enable them to be repositioned to suit the size of the equipment.

1. Place the WM-NP3 Mobile Workstation on a level floor.
2. Lock all braked castors.
3. Disconnect medical electrical equipment and remove the medical electrical equipment from shelf. If access is required, remove the medical electrical equipment from the shelf below as well.

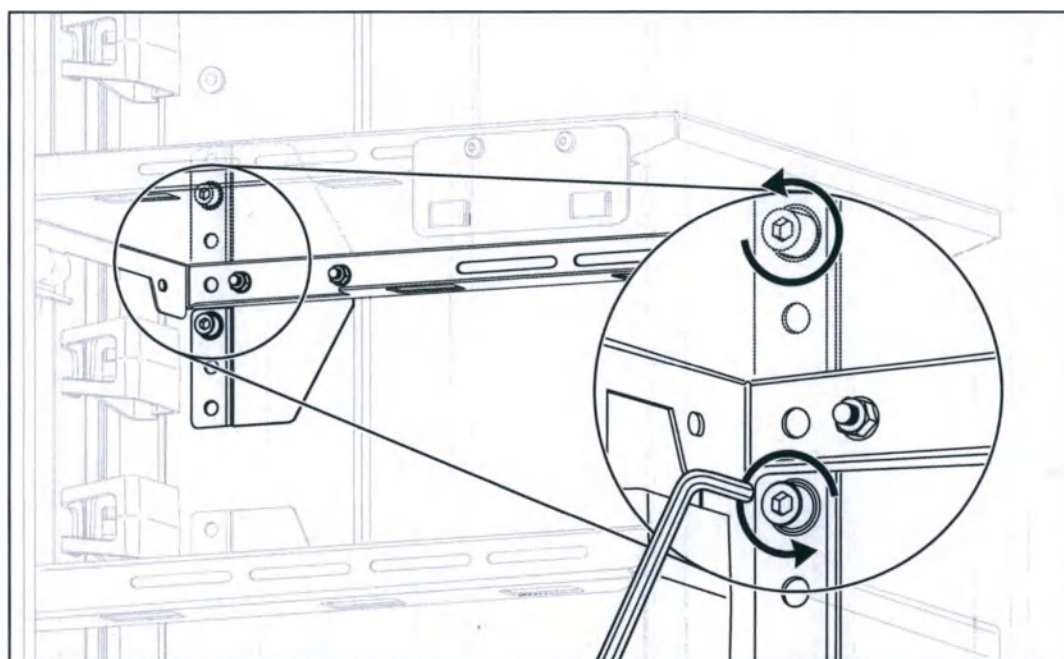4. While maintaining a firm hold on the shelf, loosen and remove the 2 M8 screws from both fixing brackets with the 6 mm hex key.
5. Move the shelf to the required height and align the holes in the bracket with the holes in the columns.
6. Insert the 2 M8 screws.
7. Check that the shelf is not tilting and tighten the 2 M8 screws (both fixing brackets) with the 6 mm hex key.

NOTICE

Risk of equipment damage

Incorrect fixing of the brackets can result in equipment damage. Always secure each bracket with 2 M8 screws.

4.1.4 MAJ-2160 Irrigation Tube Holder

The WM-NP3 Mobile Workstation is delivered with the MAJ-2160 Irrigation Tube Holder fitted to one of the intermediate shelves. It is possible to reposition this accessory in any of the slots located on the sides of the intermediate shelves.

The MAJ-2160 Irrigation Tube Holder can hold irrigation tubing with a diameter of 3.13 mm to 9.6 mm.

WARNING


Risk of infection

The MAJ-2160 Irrigation Tube Holder is not intended to hold suction tubing, as potentially infectious material from the suction tubes could contaminate the irrigation tubes and pose an infection risk to the patient.


- Do not use the irrigation tube holder to hold suction tubing.

Repositioning the MAJ-2160 Irrigation Tube Holder

1. Squeeze the clip section together from behind the intermediate shelf.

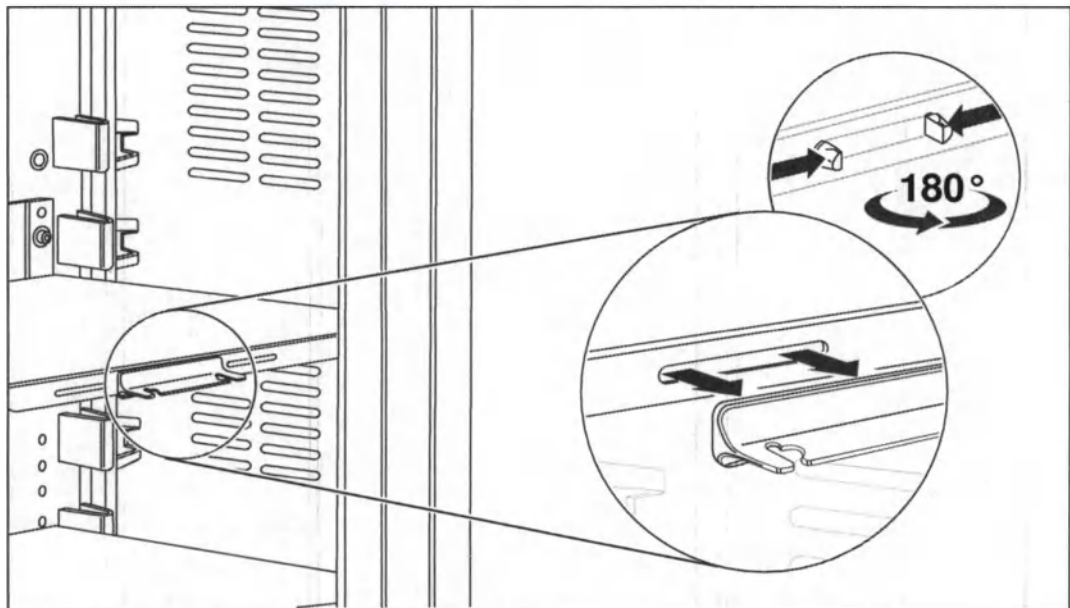
2. Pull the MAJ-2160 Irrigation Tube Holder away from the shelf to remove (view from inside).

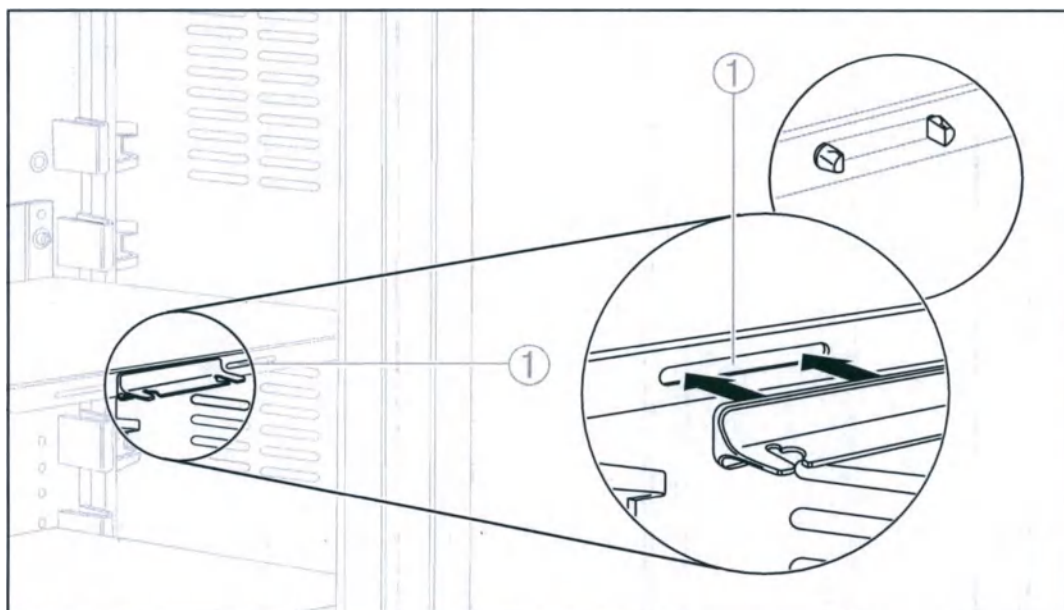3. Align the clips of the MAJ-2160 Irrigation Tube Holder with the slots (1) of the required shelf.
4. Push the MAJ-2160 Irrigation Tube Holder onto the shelf.


Using MAJ-2160 Irrigation Tube Holder

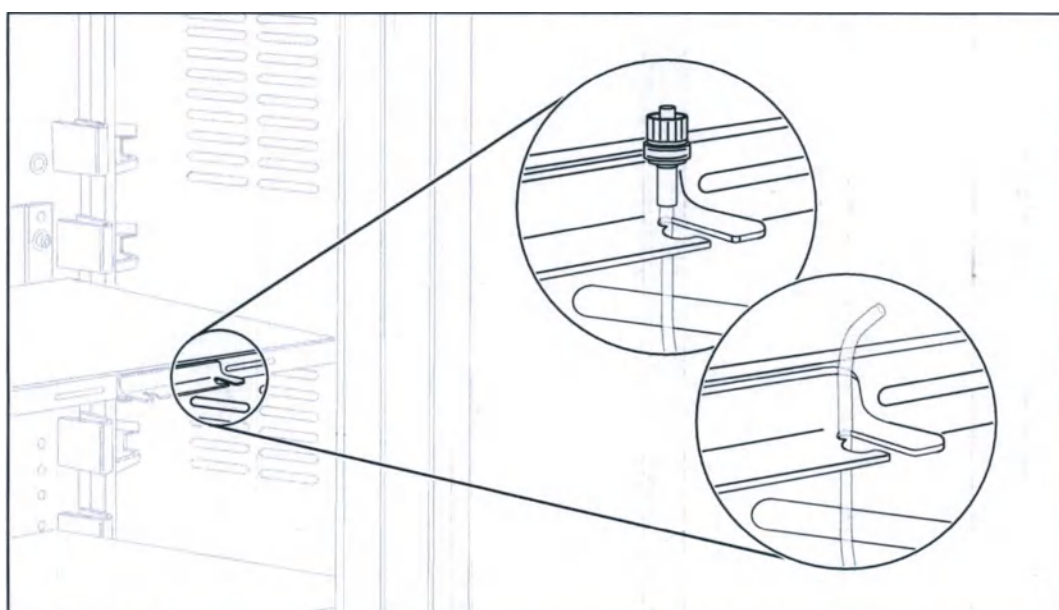- Insert the compatible irrigation tube into the irrigation tube holder by pushing it into one of the slots indicated.
- ▶ The irrigation tube is inserted.

4.1.5 Water bottle bracket

The WM-NP3 Mobile Workstation includes a pre-fitted water bottle bracket, allowing the storage of 2 water containers (not supplied).

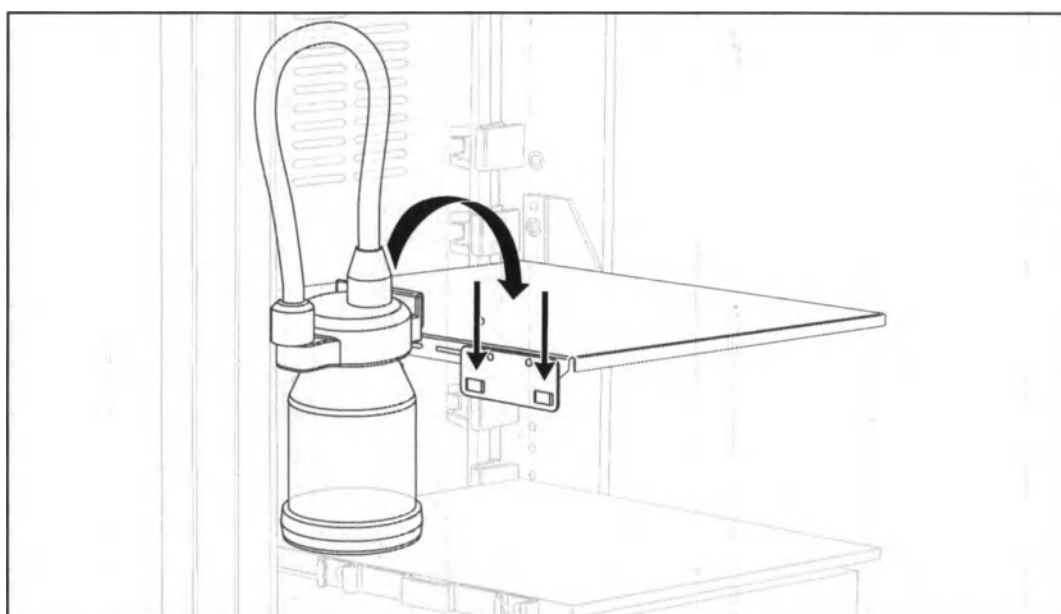The water container clips into the bracket as shown. It is possible to reposition the water bottle bracket in any of the slots located on the sides of the intermediate shelves.

4.2 Removing and attaching the rear panel

NOTICE

Risk of equipment damage

Trapping cables with the rear panel can cause cable damage. Ensure no cables are entrapped when installing and closing the rear panel.

NOTICE

Risk of equipment damage


The vent slots on the rear panel are not intended to hold equipment. Do not hang equipment on the vent slots of the rear panel.

Removing the rear panel

1. Place a hand into each of the large openings of the rear panel and pull towards you to disengage.
2. Lift the rear panel from the WM-NP3 Mobile Workstation.

Attaching the rear panel

1. Lock all braked castors.
2. To attach the rear panel, position the rear panel onto the three studs at the base of the WM-NP3 Mobile Workstation.

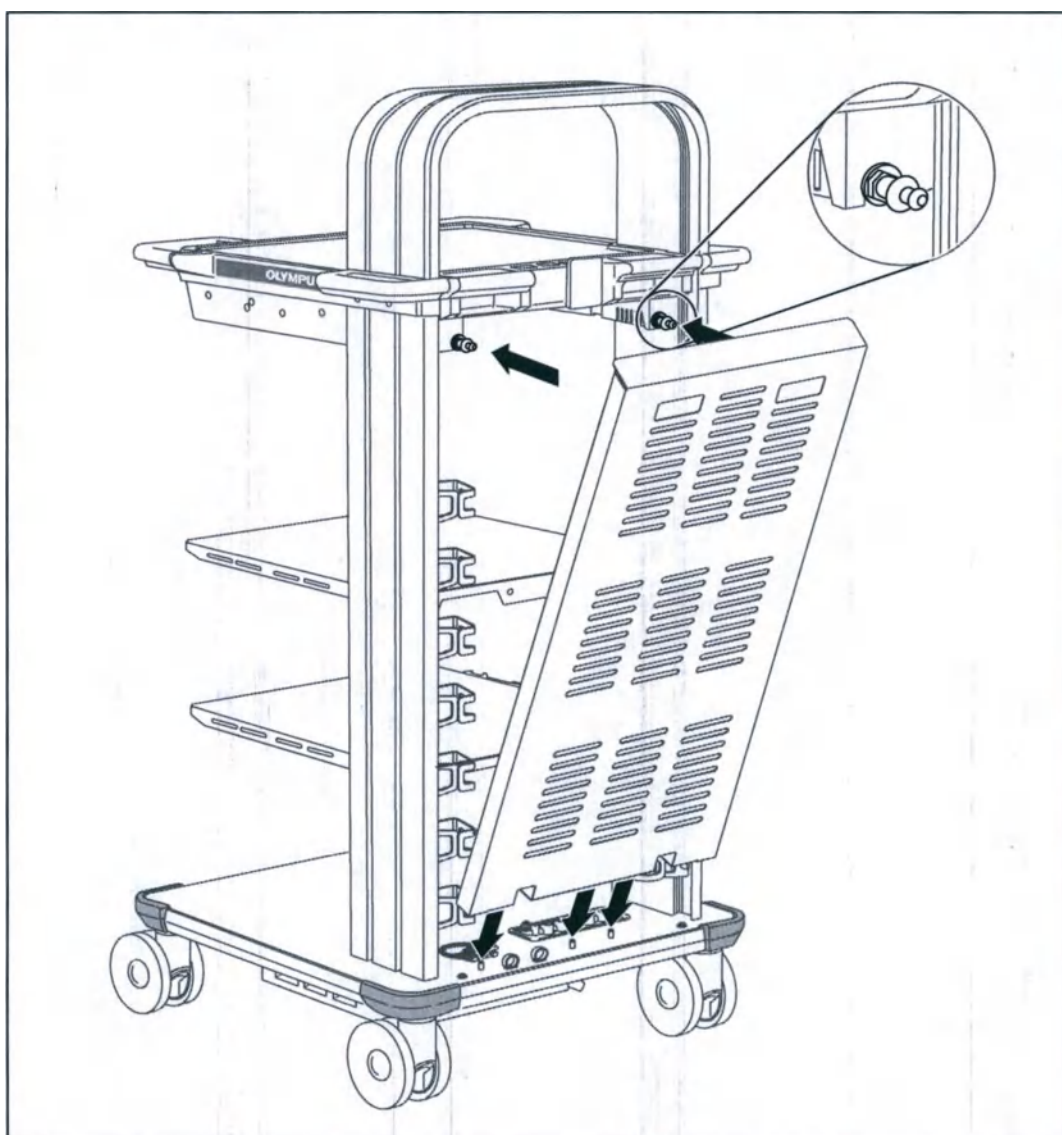
3. Push the top of the panel against the frame to engage the snaploc couplings.
4. Ensure both snaploc couplings are fully engaged.

4.3 Installing MAJ-2159 Scope Pole Kit


The MAJ-2159 Scope Pole Kit can be installed on either the left or right-hand side of the workstation. The following describes installing the MAJ-2159 Scope Pole Kit on the right hand side. Refer to “Installing MAJ-2159 Scope Pole Kit on the left-hand side” on page 41 if left-hand installation is required.

4.3.1 Installing MAJ-2159 Scope Pole Kit on the right-hand side


1. Lock all braked castors.
2. Take 4 M6 T-nuts from the fitting kit.

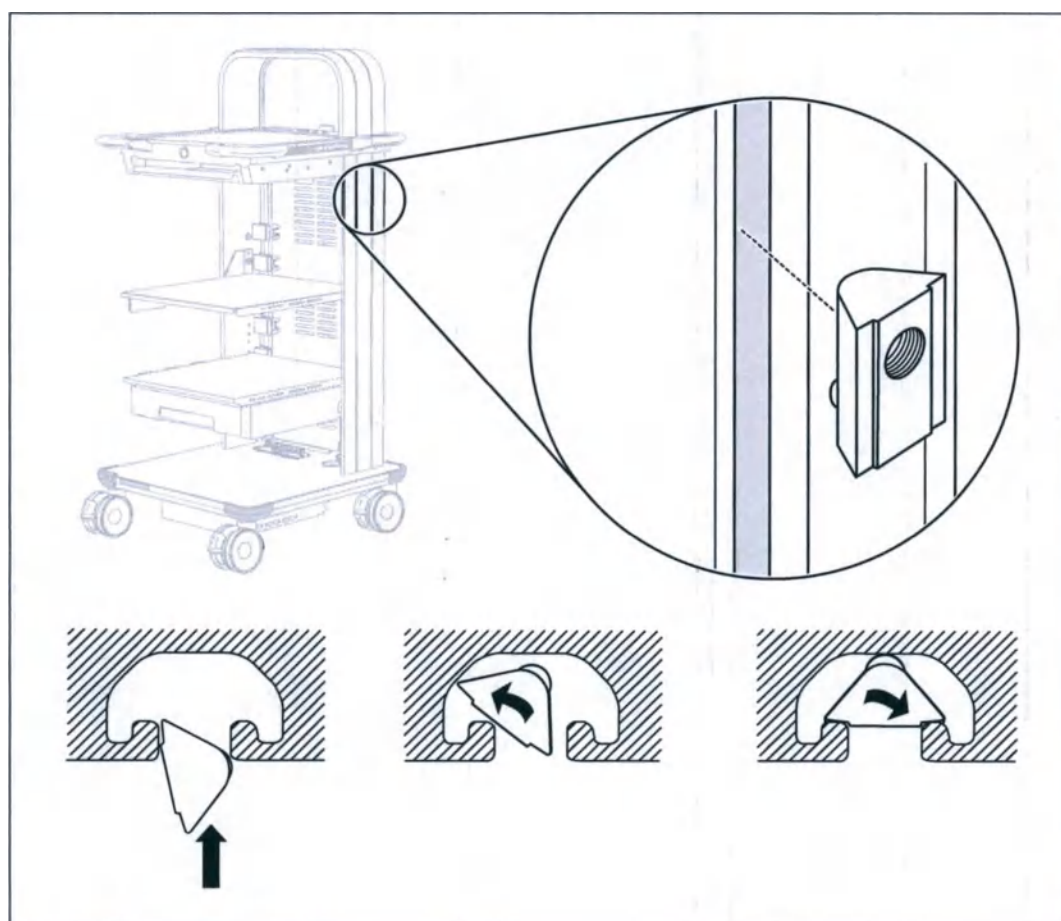
3. Insert 2 M6 T-nut, one in each slot of the WM-NP3 Mobile Workstation hoop.

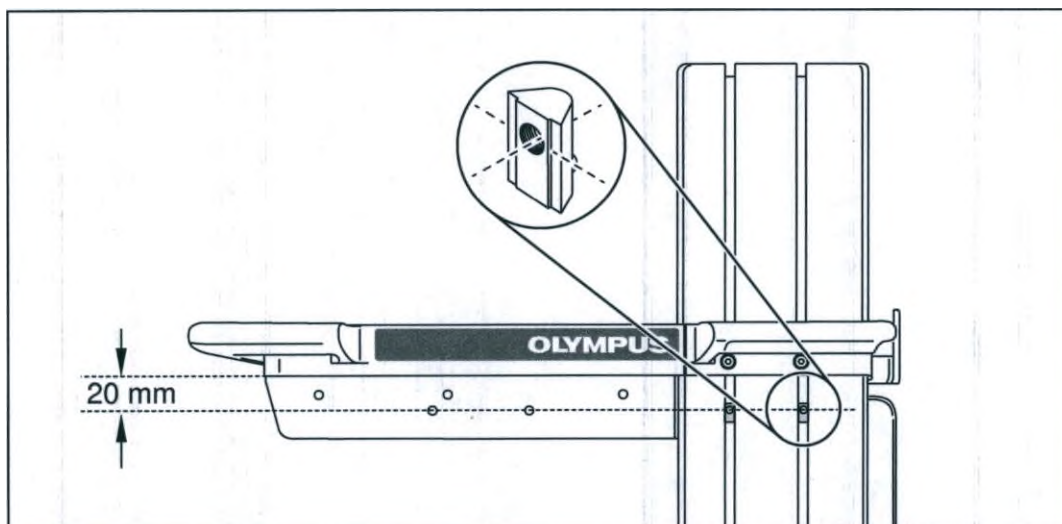
4. Position the inserted M6 T-nuts approximately 20 mm below the side handle (if fitted) or 20 mm below the top tray lower surface.

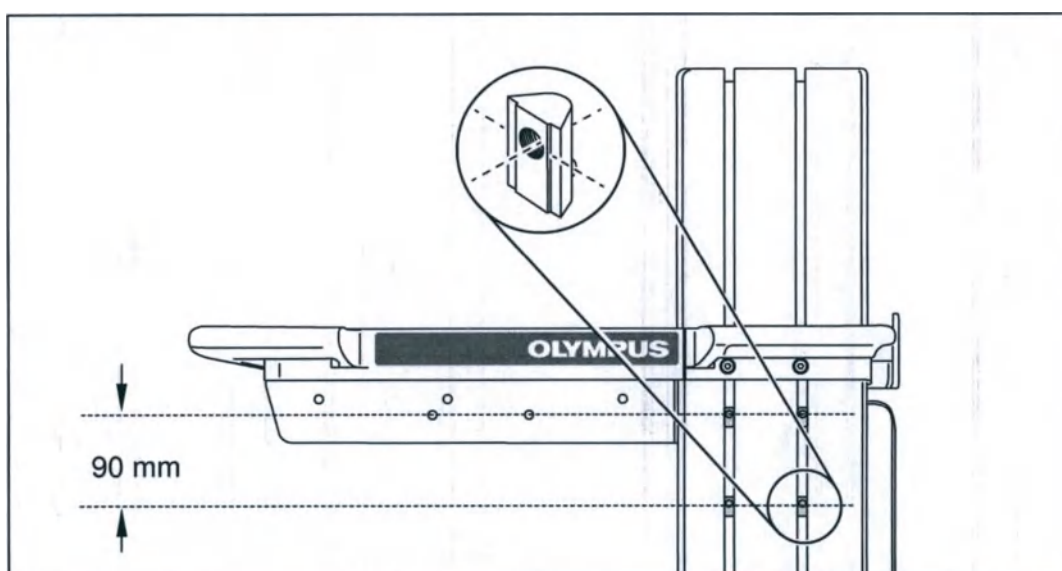5. Insert 2 more M6 T-nut, one in each slot of the WM-NP3 Mobile Workstation hoop, approximately 90 mm below the first inserted M6 T-nuts.

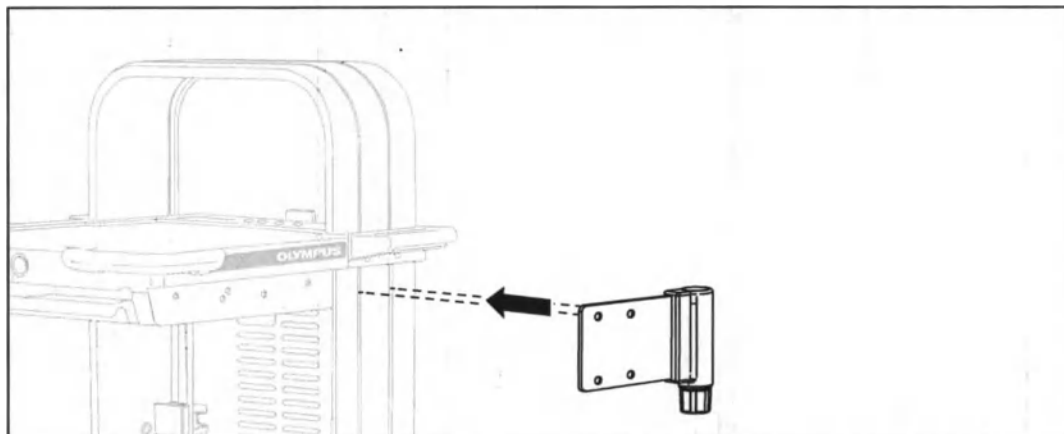6. Place the scope pole mount against the WM-NP3 Mobile Workstation hoop and align the 4 holes in the scope pole mount with the inserted M6 T-nuts.

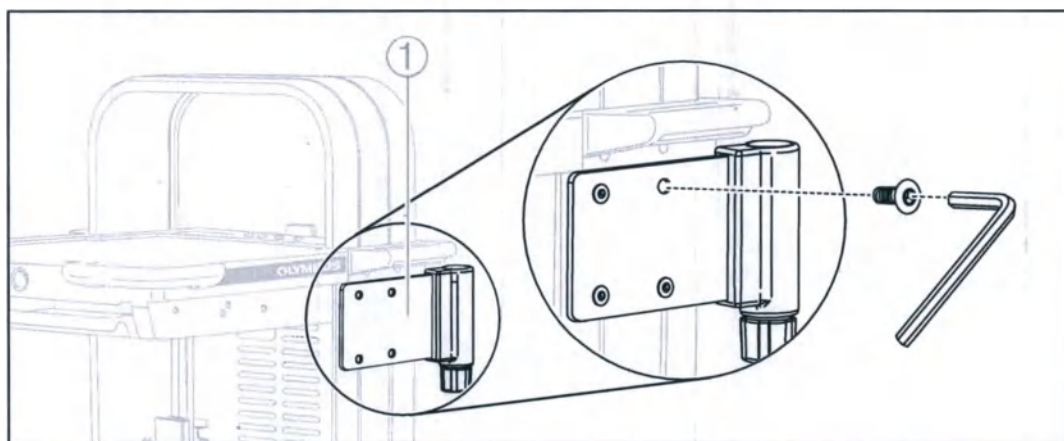- 1) Scope pole mount
7. Insert and tighten the 4 M6 countersunk screws with the 4 mm hex key to secure the scope pole mount.


8. Insert the vertical pole into the scope pole mount.

NOTICE

Risk of equipment damage

Loosening the scope pole mount lock will allow the endoscope pole to fall rapidly.

Ensure a firm grip when repositioning the endoscope pole.

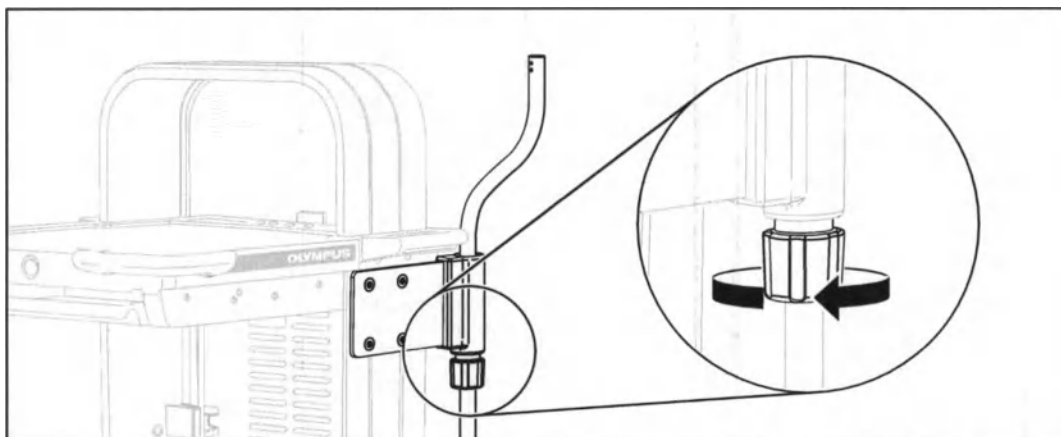
9. Ensure the vertical endoscope pole can pass through the scope pole mount by turning the scope pole mount lock.

NOTE

Do not remove the scope pole mount lock completely, unless replacement is required.

10. Move the vertical endoscope pole to the required position.

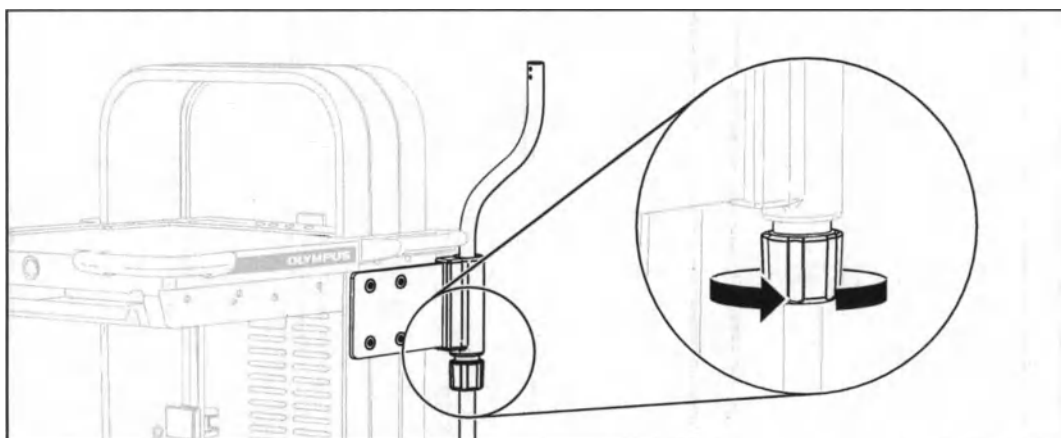
11. Lock the vertical endoscope pole by turning the scope pole mount lock.


NOTE

Turn finger tight only. Do not over-tighten.


12. Loosen and remove the pre-fitted 2 M4 screws, crinkle washers and plain washers from the vertical endoscope pole with the 2.5 mm hex key.

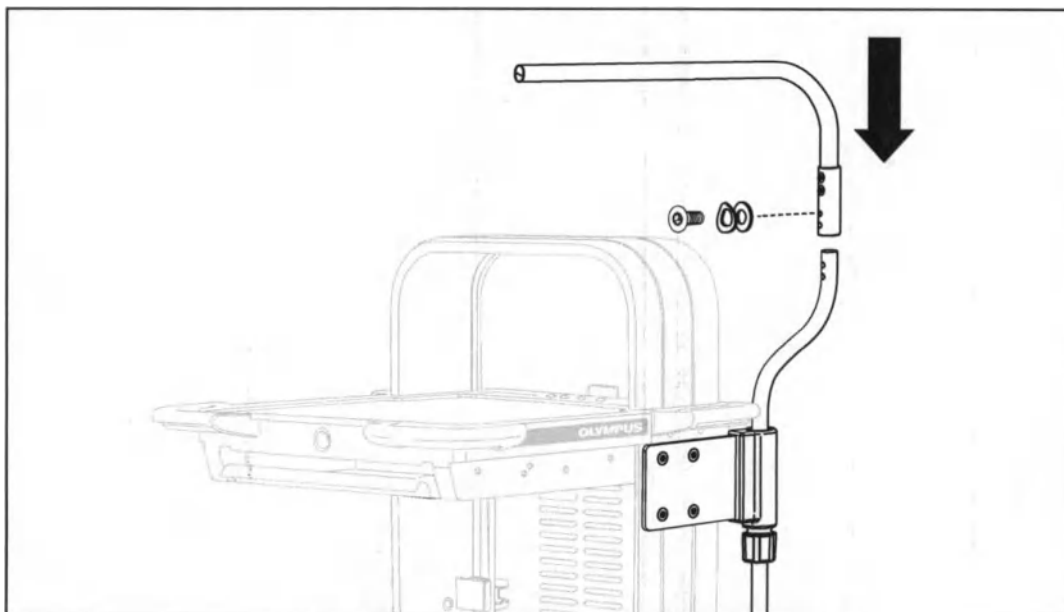
- 13. Locate the horizontal endoscope pole onto the vertical endoscope pole, align the fixing holes and tighten the 2 M4 screws, crinkle washers and plain washers with the 2.5 mm hex key.**

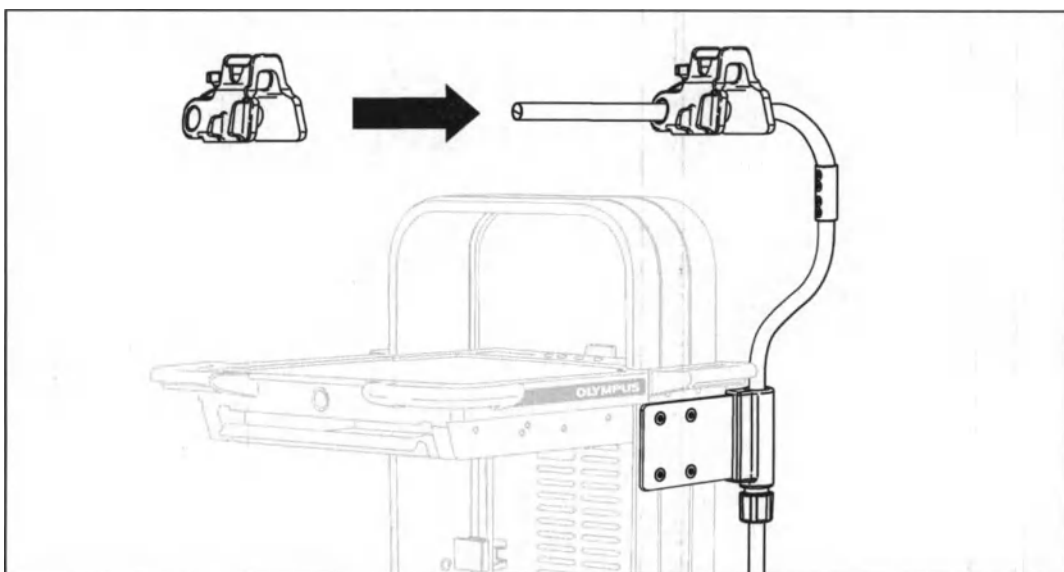
- 14. Slide the endoscope hanger onto the horizontal endoscope pole to the required position.**

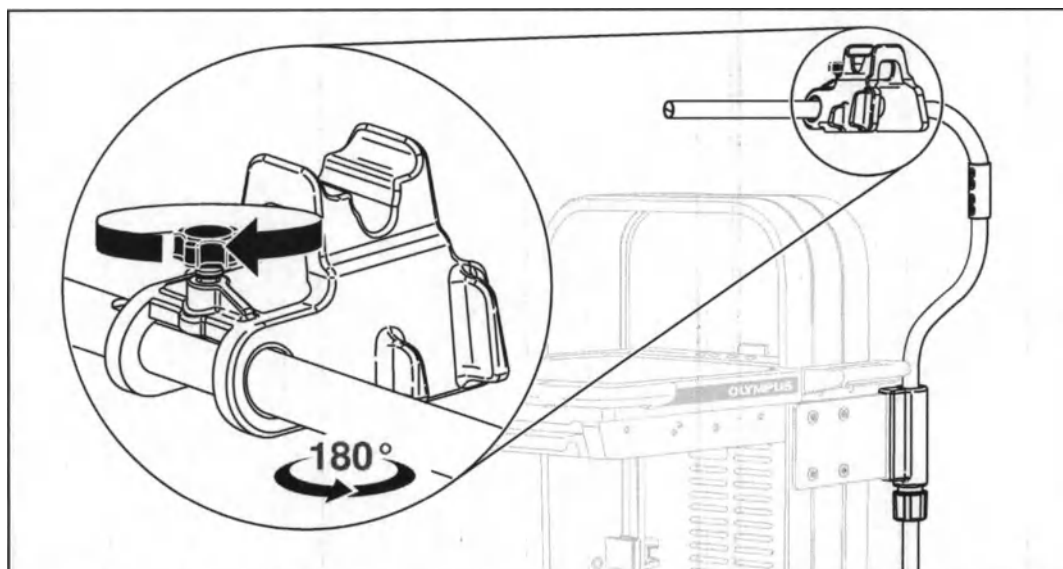
15. Tighten the handwheel to secure the endoscope hanger.
16. Repeat for second endoscope hanger.
17. To remove the endoscope hangers, loosen the handwheels and slide the endoscope hangers off the horizontal endoscope pole.

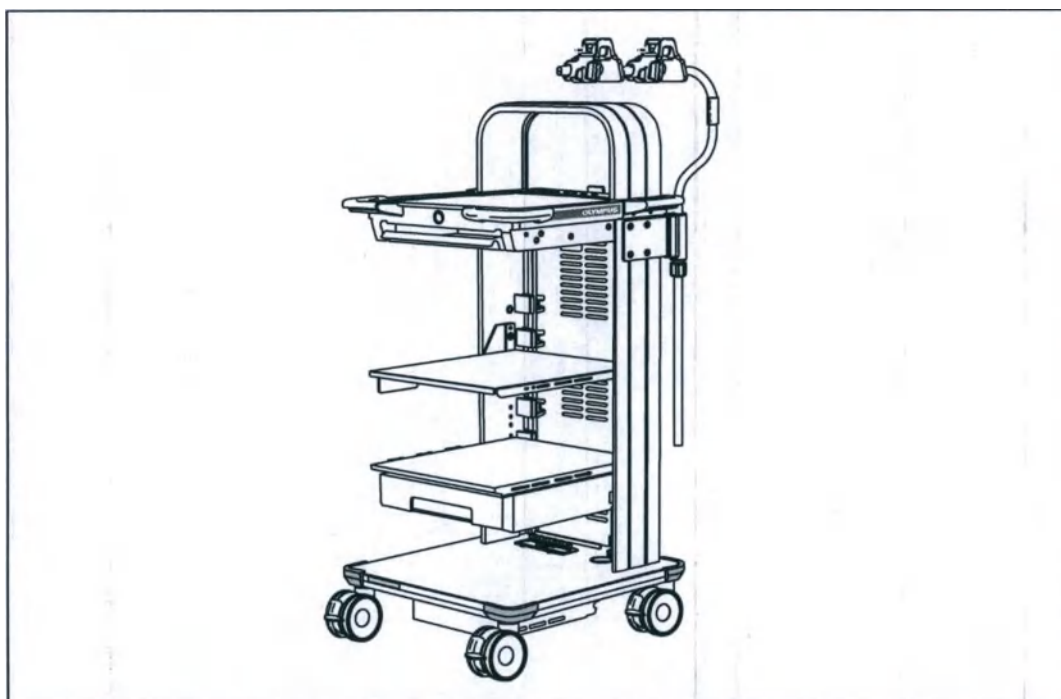18. Place the endoscope connector holder against the WM-NP3 Mobile Workstation top tray and align the 2 holes in the endoscope connector holder with the holes in the top tray.

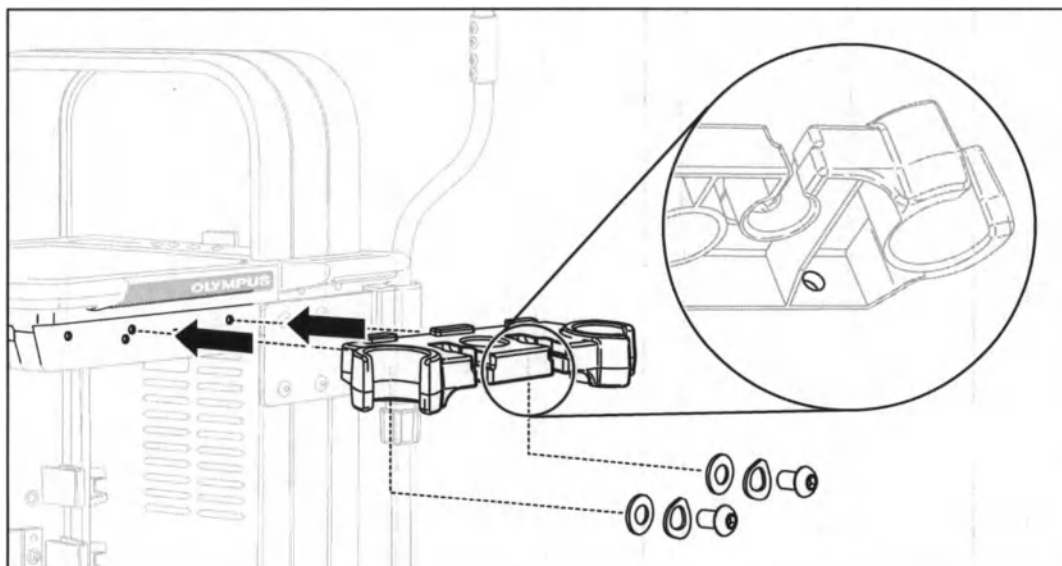- 19. Secure the endoscope connector holder with the following fastener order.**
- Position 2 M6 plain washers (1 per screw) against the endoscope connector holder, followed by 2 M6 crinkle washers.
 - Tighten 2 M6 screws (1 per screw) through the washers with the 4 mm hex key.

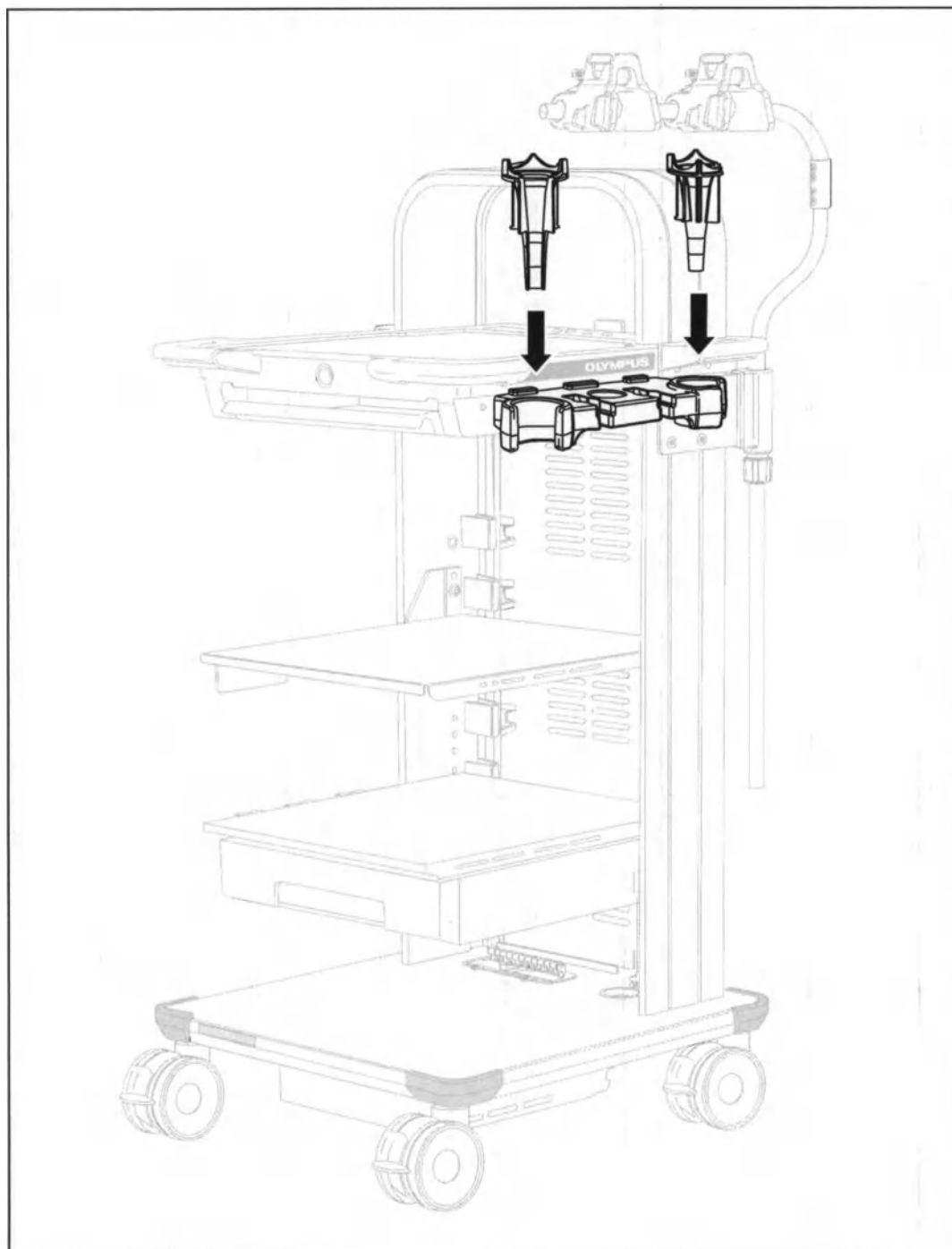20. Insert the endoscope supports into the endoscope connector holder.


NOTE

The endoscope connector holder also provides an area to fit a white balance cap (not supplied).

4.3.2 Installing MAJ-2159 Scope Pole Kit on the left-hand side

For the left-hand installation of MAJ-2159 Scope Pole Kit to the workstation, the scope pole mount installation needs to be changed. If the scope pole mount is already installed on the right-hand side, the scope pole mount needs to be removed from the workstation first.


Reversing the scope pole mount for left-hand installation

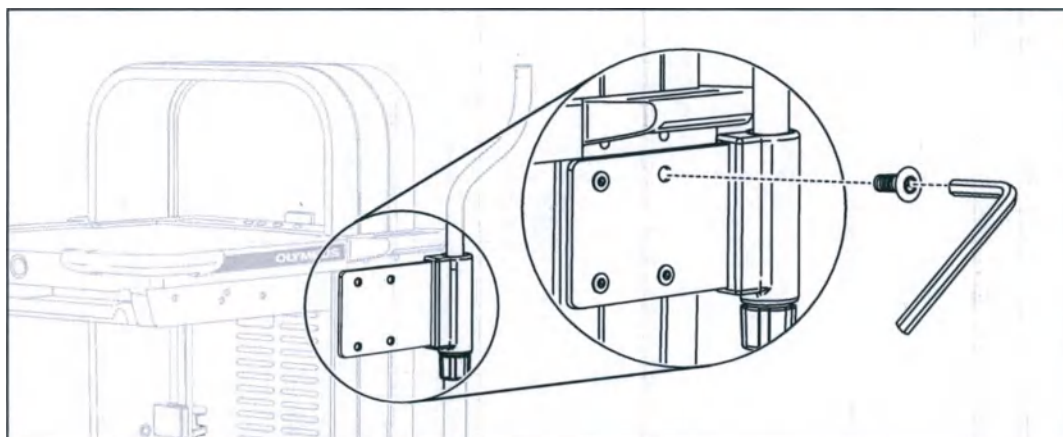
1. Loosen the 4 M6 countersunk screws with the 4 mm hex key.


2. Use a small flat-blade screwdriver and remove the 4 M6 T-nuts from the slots in the WM-NP3 Mobile Workstation hoop.

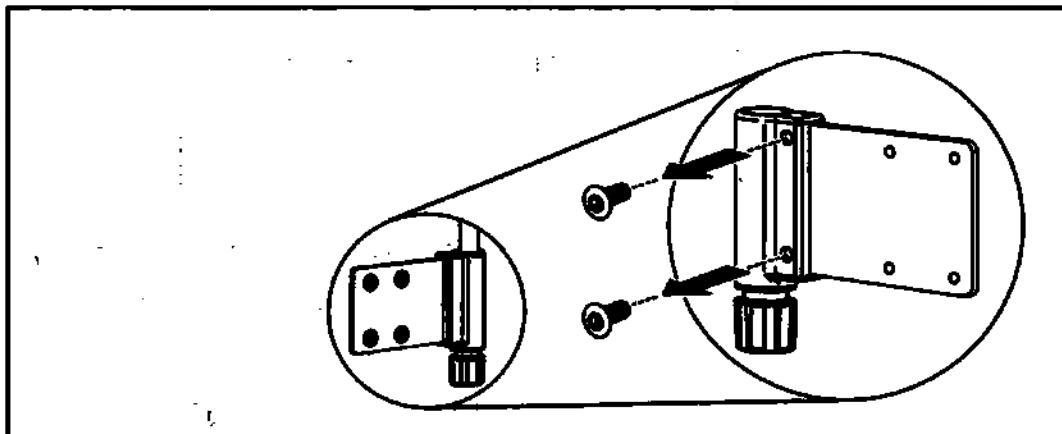
3. Loosen the 2 M8 countersunk screws from the scope pole mount with the 5 mm hex key.

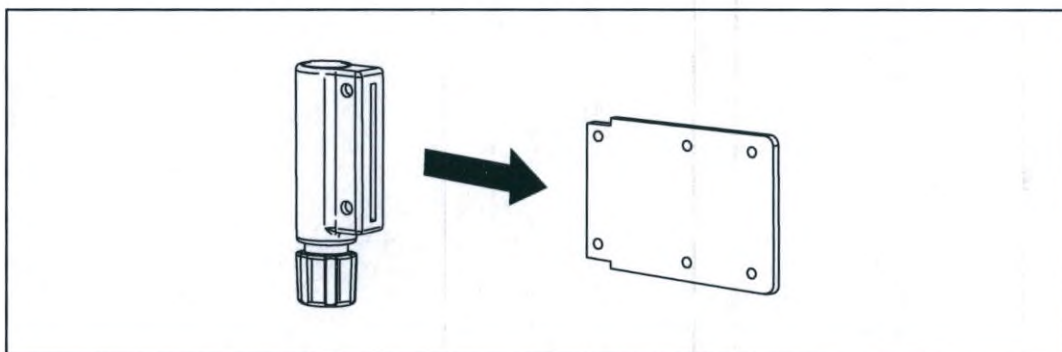4. Disassemble the scope pole mount.

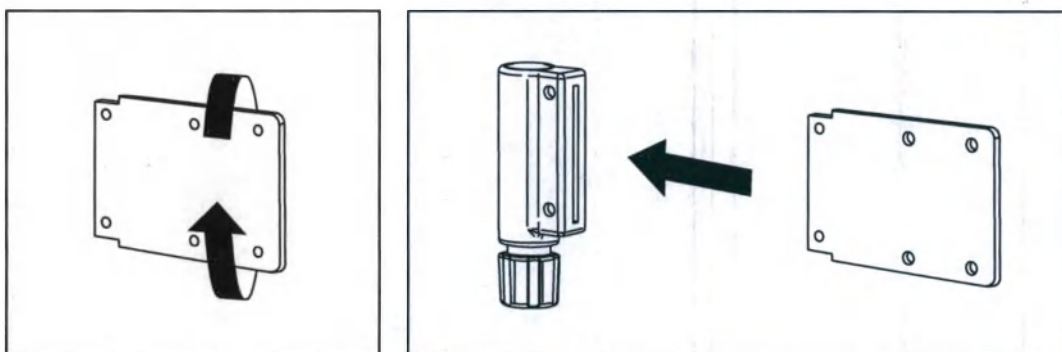5. Turn the plate over and re-insert it.

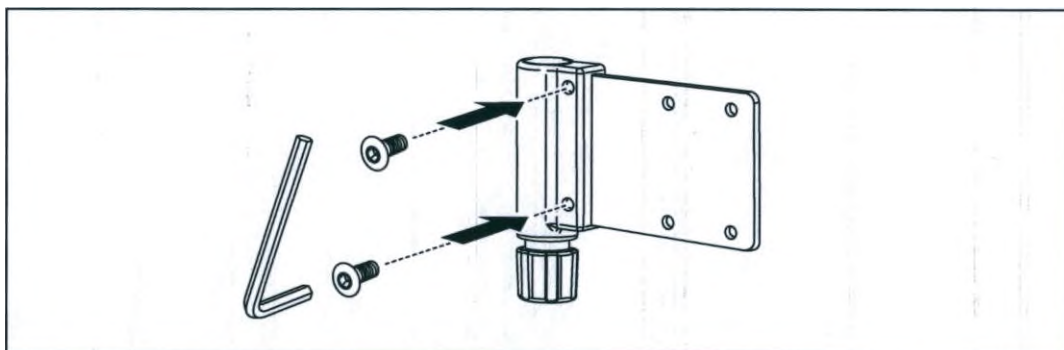
6. Re-tighten the 2 M8 countersunk screws with the 5 mm hex key.


The scope pole mount can now be installed on the left-hand side of the WM-NP3 Mobile Workstation. Refer to the instructions in "Installing MAJ-2159 Scope Pole Kit on the right-hand side" on page 32.

4.4 Installing keyboards

1. Fully extend the keyboard tray.

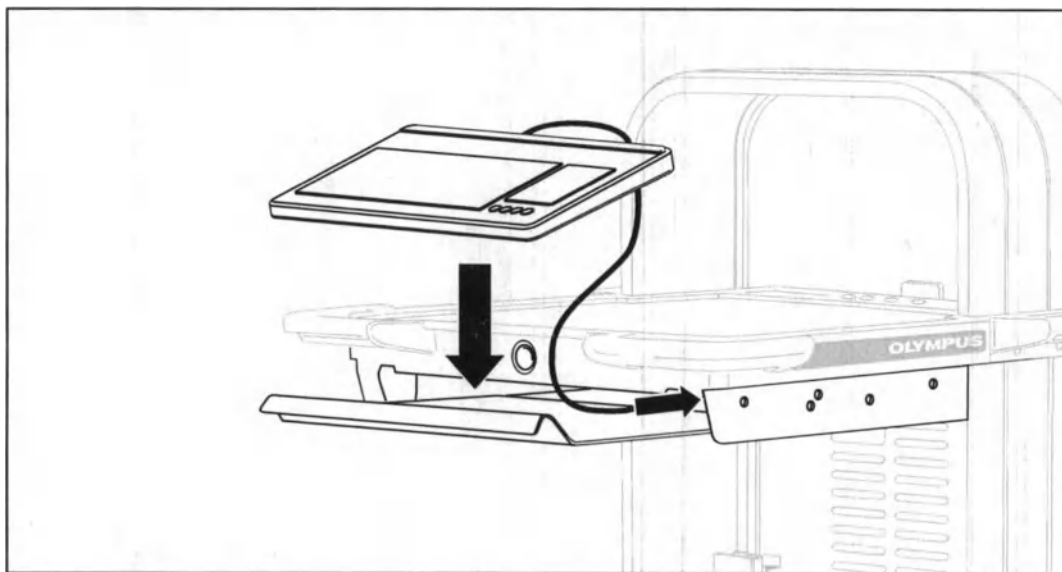2. Place the keyboard on the keyboard tray.
3. Remove the rear panel if necessary, refer to "Removing the rear panel" on page 31.

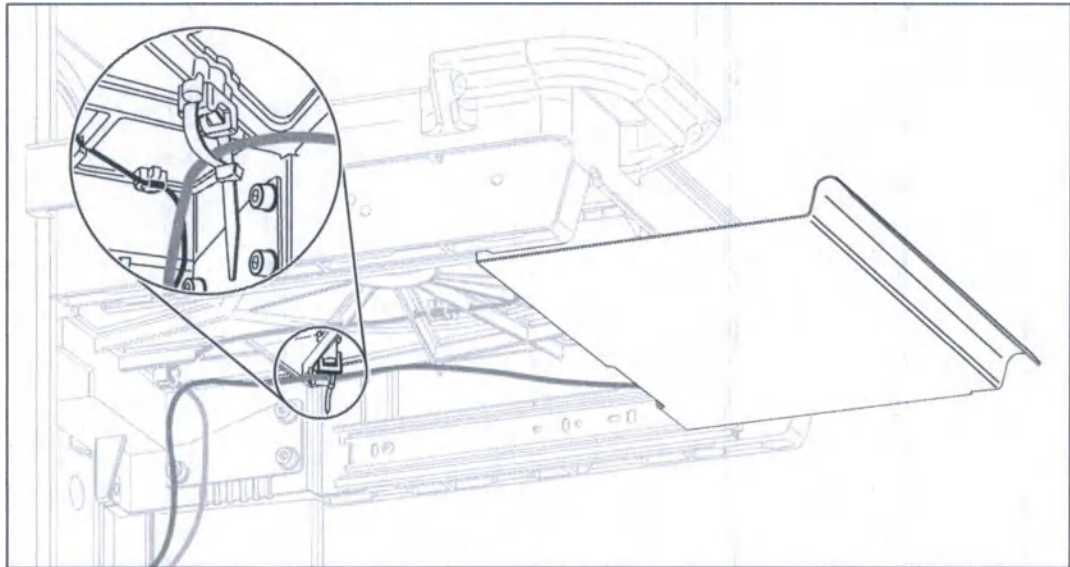
4. Route and secure the keyboard cable with a cable tie in the cable clip on the underside of the top tray.
5. **Connect the keyboard cable to the medical electrical equipment.**
6. **Re-attach the rear panel.**

4.5 Installing medical electrical equipment

Ensure the fully installed system meets the requirements described in “Safety standards for medical electric systems” on page 71. For advice, contact your Olympus representative or qualified service personnel.

4.5.1 View of the separation transformer

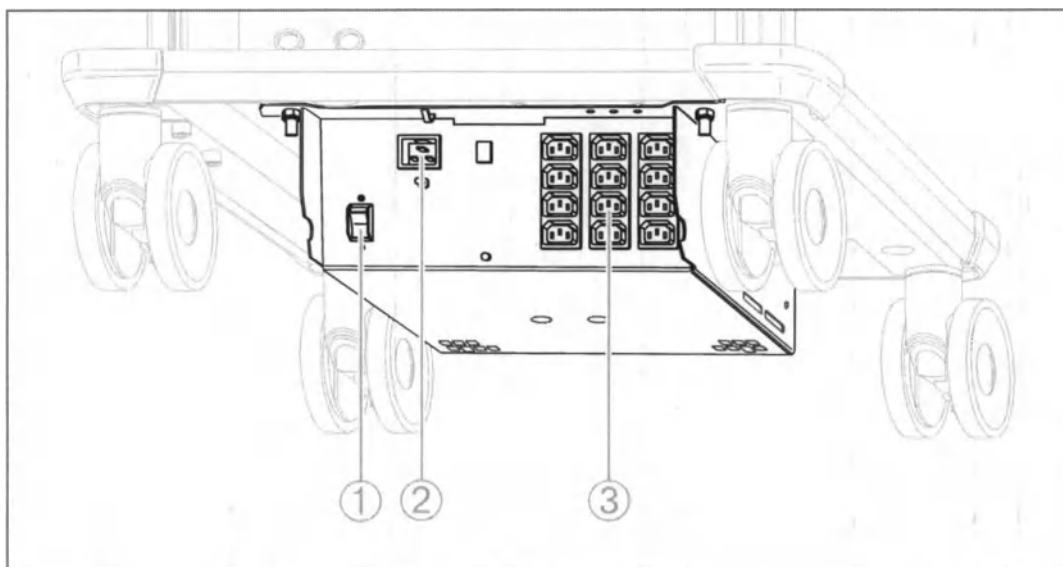- 1) Power switch / circuit breaker
- 2) Input socket for power cable
- 3) Output-sockets for medical electrical equipment

4.5.2 Installing medical electrical equipment

CAUTION

Risk of injury

Incorrect work practices when installing medical electrical equipment on the WM-NP3 Mobile Workstation can result in equipment damage. Observe the following.


- Place medical electrical equipment on the WM-NP3 Mobile Workstation by 2 persons working together.
- Do not exceed the load capacity detailed in “Technical data” on page 66.
- Only use the supplied power cable of the medical electrical equipment.

NOTICE

Risk of equipment damage


Ensure the medical electrical equipment connected to the separation transformer does not exceed the separation transformer's rating.

Installing medical electrical equipment

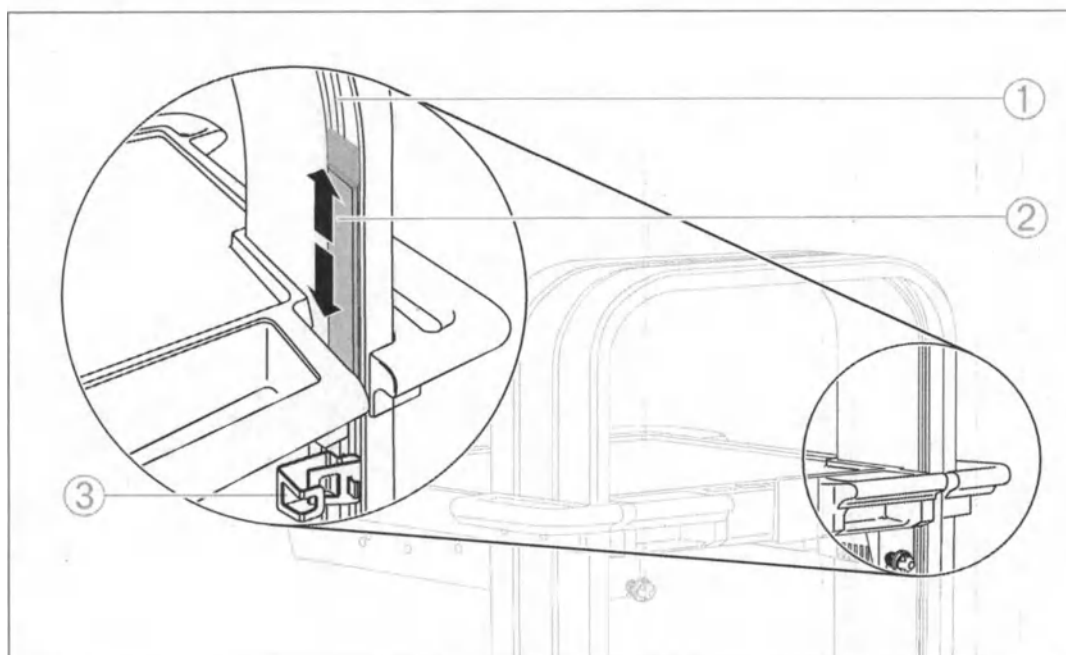- 1) Slot
- 2) Cable management strip
- 3) Cable management unit
1. Check that the medical electrical equipment is switched off.
2. Check that the connector of the separation Transformer power switch cable is secure.
3. Place medical electrical equipment on the shelves, referring to "Maximum size of equipment for shelves" on page 24.
4. Make sure medical electrical equipment is placed centrally on the shelves.

5. Use the cable management features of the workstation hoop to store and organize cables:
 - For cables above the top tray (e.g. for the LCD monitor arms), store the cables behind the cable strip that can be pushed around the workstation hoop. Pull the cable strip forward to access the slot and push it back in when the cables are positioned.
 - For cables below the top tray, lead the cables through the cable management units in the lower areas of the workstation.
6. Connect the power cable of the medical electrical equipment to the separation transformer. Ensure all connectors are fully inserted.
7. Secure power cables with cable ties from the fitting kit, without loosening the power cables from the separation transformer output-sockets.

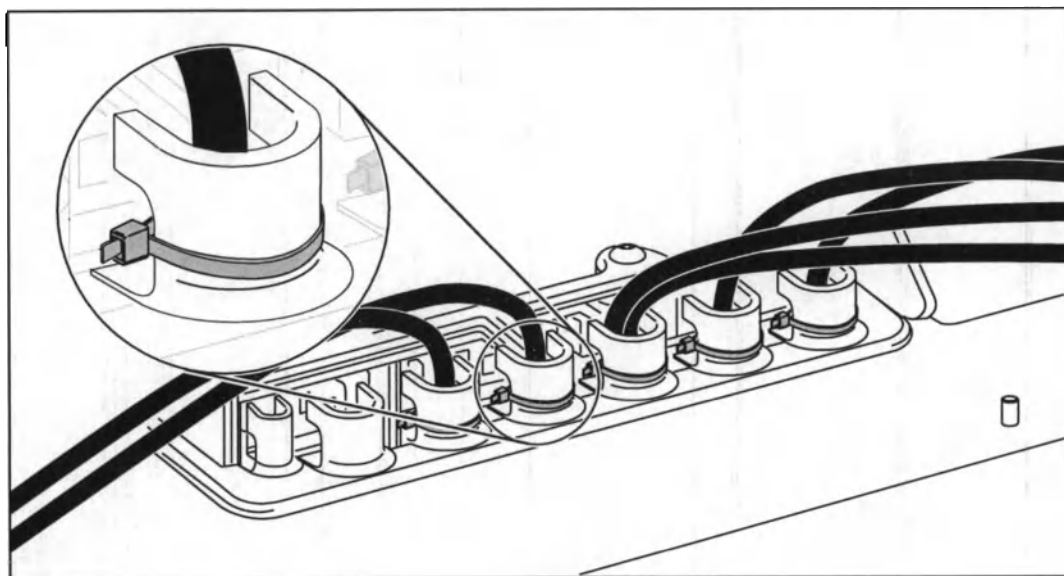NOTE

Do not pull the cables taut. This could result in excessive strain on the cable and/or connector interface.

4.5.3 Connecting medical electrical equipment to electricity

**CAUTION****Risk of injury**

If a failure of the power supply would subject the patient to an unacceptable risk, supply the workstation from an uninterruptible power supply (UPS) or a supply with a back-up generator.


NOTICE**Risk of equipment damage**

Taut cables can result in equipment damage. Ensure that all power cables have sufficient slack so that they cannot be pulled taut when moving the WM-NP3 Mobile Workstation around the working area. Relocate if necessary.

**NOTE**

Operating the power button several times in rapid succession can cause the mains power fuse to operate. If this should occur, reset the circuit breaker(s) by pressing in the power switch / circuit breaker button as shown in section “4.5.1 View of the separation transformer” on page 45. See also section “12 Troubleshooting” on page 99.

1. Connect the WM-NP3 Mobile Workstation power cable and light source power cable (where separate) to a properly grounded mains electricity. The power outlets should be easily accessible when the mobile workstation is in its normal operating position.

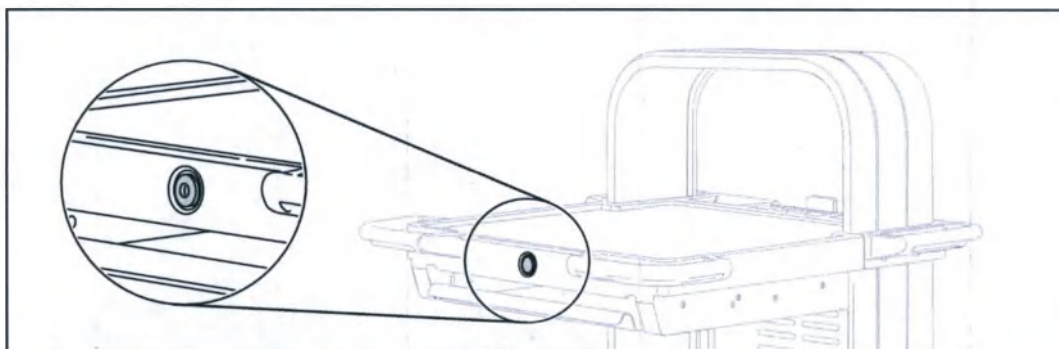- The power button on the WM-NP3 Mobile Workstation illuminates white.

Before use

- ▶ The LED on the separation transformer illuminates green.
2. Switch on the WM-NP3 Mobile Workstation with the power button.
 - ▶ The power button illuminates green.
3. Switch on the connected medical electrical equipment.
4. Check that all installed medical electrical equipment operates correctly.
5. Attach the rear panel, taking care not to trap any cables.

5 Maneuvering

This chapter provides information on maneuvering the product.

5.1 Specific safety notes for maneuvering

WARNING

Risk of explosion

None of the medical electric equipment with which the WM-NP3 Mobile Workstation is intended to be used is 'AP' or 'APG' rated.

- Do not install the WM-NP3 Mobile Workstation within the zone of risk of flammable gases.

CAUTION

Risk of injury

CO₂ cylinders weighing more than 29 kg can cause instability during maneuvering, resulting in injury. Always remove CO₂ cylinders weighing more than 29 kg before maneuvering.

NOTICE

Risk of equipment damage

Incorrect work practices when maneuvering the WM-NP3 Mobile Workstation can result in equipment damage. Observe the following.

- Use 2 persons to maneuver the WM-NP3 Mobile Workstation in between rooms.
- Take particular care when maneuvering the WM-NP3 Mobile Workstation up or down slopes.


5.2 Preparing for maneuvering

1. Switch off the WM-NP3 Mobile Workstation.
2. Disconnect the power cable from the mains electricity.
3. Hang the power cable on the cable hook on the back of the top tray.


4. Remove CO₂ cylinders weighing more than 29 kg and transport the cylinders separately. If the cylinder weight is unknown, confirm the weight with the hospital's gas cylinder supplier.
5. Check that the installed medical electrical equipment is stable.
6. Put the accessories into their respective storage positions.
7. Place the LCD Monitor Arm (if fitted) in the retained position. If the LCD Monitor Arm blocks your view while maneuvering, rotate the LCD Monitor Arm (still in the retained position) towards the front of the WM-NP3 Workstation.

5.3 Maneuvering the WM-NP3 Mobile Workstation

1. Unlock all braked castors.
2. Move the WM-NP3 Mobile Workstation only by gripping the handles on the front, sides and back of the top tray. The WM-NP3 Mobile Workstation can also be maneuvered with the side handles on the hoop (if fitted).

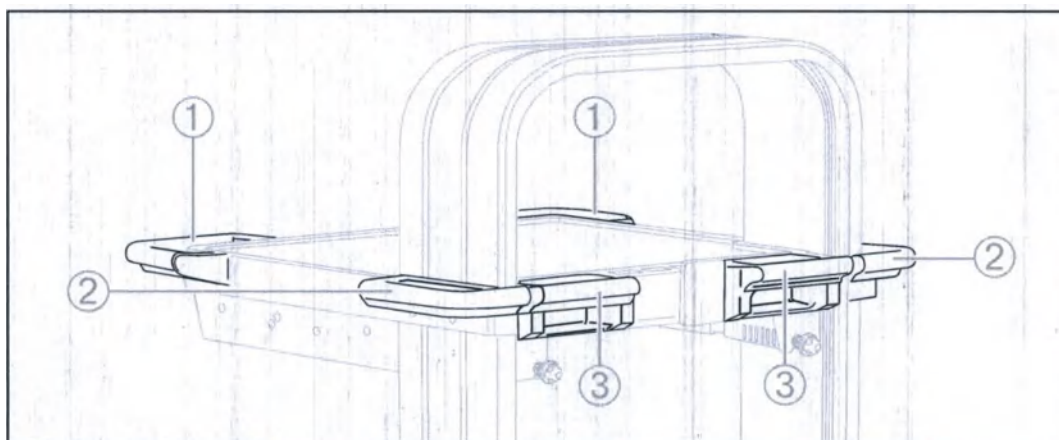- 1) Front handle
- 2) Side handle
- 3) Back handle

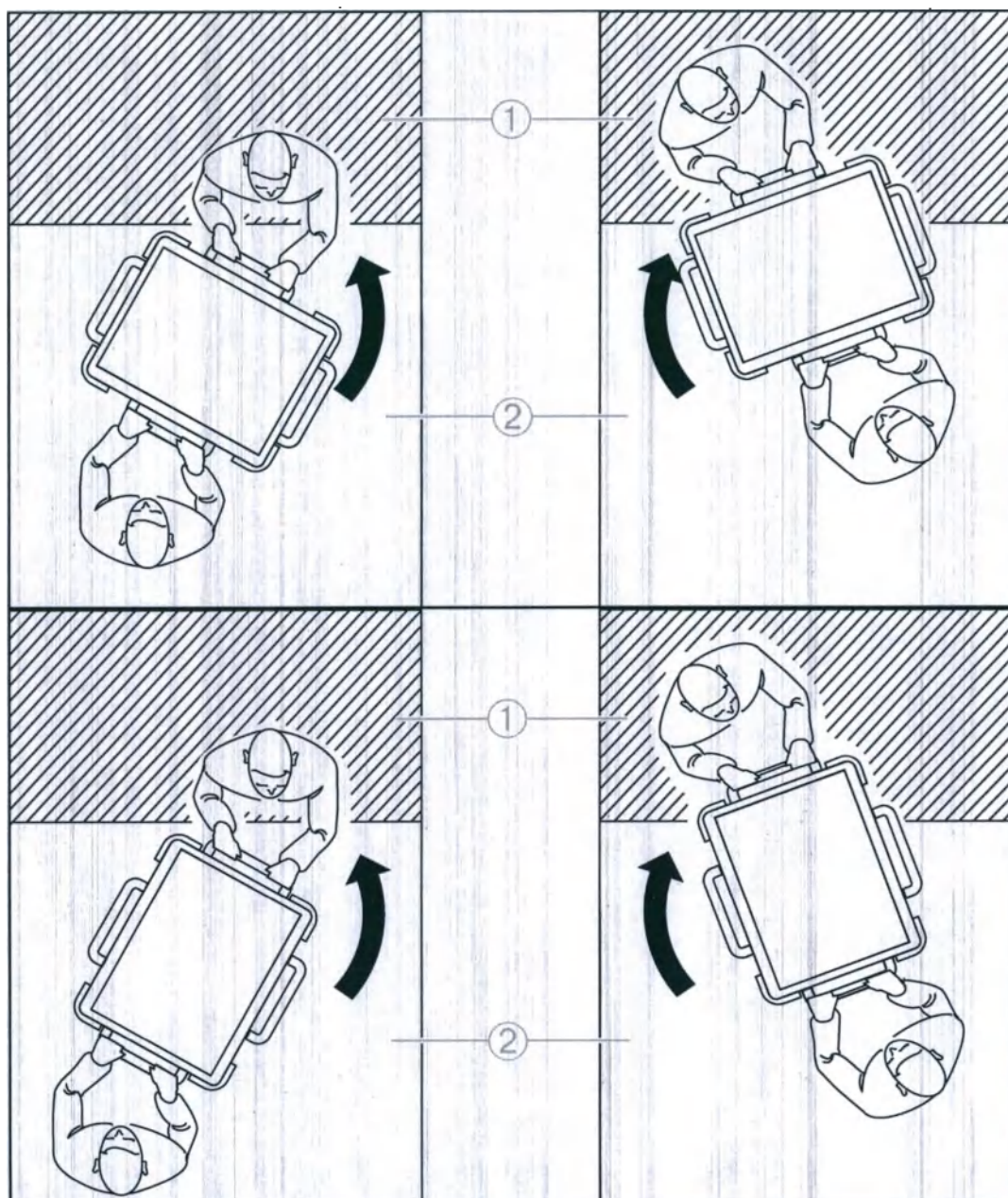- 1) Higher level
- 2) Lower level

- 3. When crossing cables or doorsteps, move the product slowly so that the castors pass over the obstacle one at a time.
- 4. After maneuvering, place the WM-NP3 Mobile Workstation on a level floor.
- 5. Lock all braked castors.

6 Use

This chapter provides information on using the product.

6.1 Specific safety notes for use

WARNING

Risk of explosion

None of the medical electric equipment with which the WM-NP3 Mobile Workstation is intended to be used is 'AP' or 'APG' rated.

- Do not use the WM-NP3 Mobile Workstation within the zone of risk of flammable gases.

WARNING

Risk of electric shock

When applied directly to the patient's heart, the combination of medical electrical equipment used on the WM-NP3 Mobile Workstation can result in patient shock.

- Do not apply the combination of medical electrical equipment of the WM-NP3 Mobile Workstation directly to the patient's heart.

WARNING

Risk of electric shock

Do not modify this equipment without authorisation of the manufacturer.

**WARNING****Risk of infection**

Incorrect work practices when handling the WM-NP3 Mobile Workstation can increase the risk of infection.

- Do not hang contaminated endoscopes on the endoscope pole.

**CAUTION****Risk of injury to the patient and user**

Incorrect work practices when handling installed medical electrical equipment can result in injury to the patient.

- Check that all installed medical electrical equipment operates correctly before using the system 'on patient'.

NOTICE**Risk of equipment damage**

Incorrect work practices when handling liquids around the WM-NP3 Mobile Workstation can result in equipment damage.

Observe the following.

- Keep open containers of liquid away from the WM-NP3 Mobile Workstation.
- Do not use the WM-NP3 Mobile Workstation if liquid has entered the system. Immediately isolate the WM-NP3 Mobile Workstation from the mains supply by disconnecting the mains power plug (refer to section "10.1 Power cable" on page 66) and contact an Olympus representative.


NOTICE**Risk of equipment damage**

Incorrect work practices when operating the WM-NP3 Mobile Workstation can result in equipment damage. Observe the following.


- Do not lean or stand on any part of the WM-NP3 Mobile Workstation or accessory.
- Do not place any items on top of a monitor (if fitted).
- Do not exceed the maximum storage.

6.2 Using the WM-NP3 Mobile Workstation

1. Place the WM-NP3 Mobile Workstation on a level floor.
2. Ensure the WM-NP3 is positioned at least 20 cm away from any walls.
3. Lock all braked castors.
4. Connect the power cable of the WM-NP3 Mobile Workstation to the mains electricity.
 - ▶ The power button illuminates white.
5. Switch on the WM-NP3 Mobile Workstation.
 - ▶ The power button illuminates green.

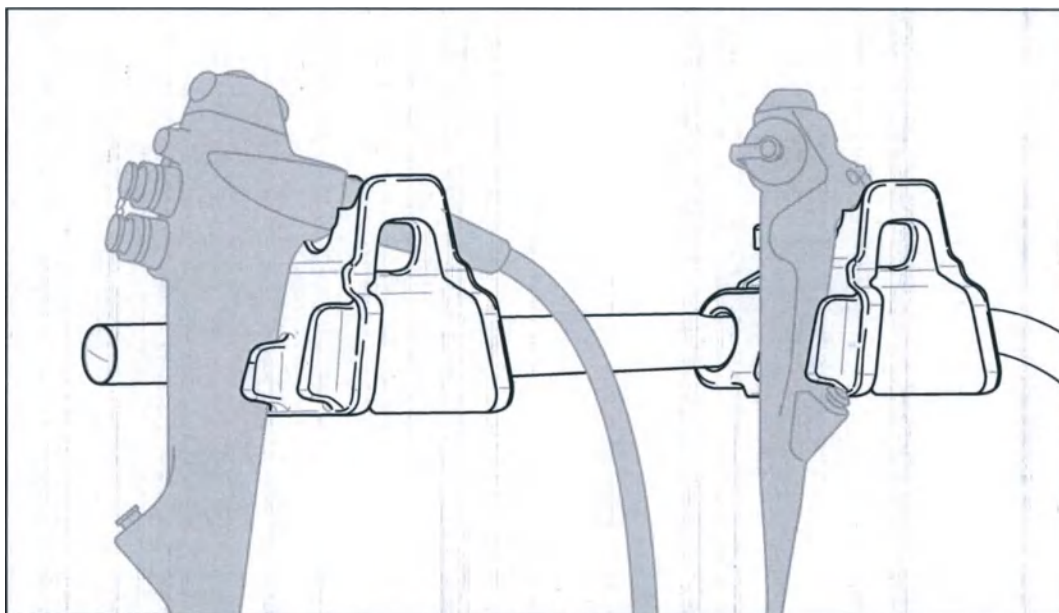6. Place the endoscopes on the endoscope hangers.

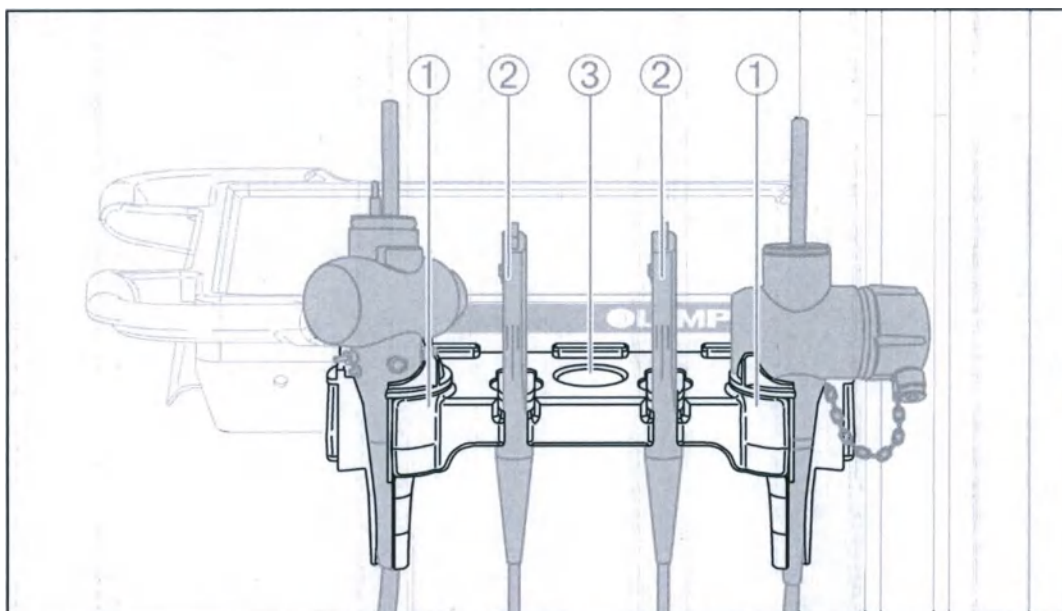7. The endoscope connector holder can store the following.
 - Endoscope connectors in the endoscope supports (1).
 - Card edge connectors in the respective cutout (2).
 - A white balance cap in the respective cutout (3).
8. For use of equipment and accessories, refer to the respective instructions for use.
9. After use, remove endoscopes and any water container from the WM-NP3 Mobile Workstation.
10. Switch off the WM-NP3 Mobile Workstation after use.
 - The power button illuminates white.
11. If the WM-NP3 Mobile Workstation will not be used for a significant period of time or if it will start the reprocessing cycle, disconnect the power cable and hang the power cable on the cable hook on the back of the top tray.

**NOTE**

Even in 'stand-by' condition, the separation transformer consumes a small amount of power, less than 2% of its rating.

7 Reprocessing

Brand-new products must be reprocessed with a complete reprocessing cycle. After each use, immediately perform the reprocessing procedures as described in the following sections. If cleaning is delayed, debris encrustation can become a source of infection.

7.1 Safety notes for reprocessing

WARNING

Risk of infection

Potentially infectious material from the patient can pose an infection risk. Reprocess the product according to the instructions in this chapter if the equipment becomes soiled with blood or other potentially infectious materials.

NOTICE

Risk of equipment damage

Wiping the surface of the WM-NP3 Mobile Workstation with abrasive cloths or abrasive detergents can cause material deterioration and surface stressing, causing paint to lift and dirt traps form. Only use the cleaning and disinfection methods as described in this chapter.

7.2 Using reprocessing agents

The reprocessing agent must be recommended by the manufacturer for cleaning and disinfection of medical devices and for the material to be reprocessed.

The reprocessing agent must be approved and validated by the FDA, EPA or a comparable institution for its efficiency.

Using reprocessing agents

- Refer to national and local guidelines regarding the approval of using alcohol as disinfectant.

- Select contact times and concentrations of the solution according to the instructions given by the manufacturer of the disinfectant.

The WM-NP3 Mobile Workstation and its ancillary components have been chemically verified to withstand the following cleaning agents.

- 2% aqueous neutral detergent
- 70% ethanol
- Benzalkonium chloride (alkyldimethylbenzylammonium chloride)
- Sodium hypochlorite NaClO (see dilution table below)

Sodium hypochlorite NaClO (CAS#7681-52-9) active ingredient	Dilution factor	
	# Parts water	# Parts NaClO
5.25%	9	1
6.15%	10	1
8.25%	14	1

Confirm the use of cleaning chemicals not listed with your Olympus service center prior to use.

Please note that response time will depend on the type of chemical being considered for use.

7.3 Cleaning and disinfection

7.3.1 Safety notes for cleaning and disinfection

WARNING

Risk of infection

Incorrect work practices when cleaning and disinfecting the WM-NP3 Mobile Workstation can increase the risk of infection. Observe the following.

- Use personal protective equipment to guard against dangerous chemicals and potentially infectious material.
- Wear the following protective equipment.
 - Eye wear
 - Face mask
 - Moisture-resistant clothing
 - Chemical-resistant gloves that fit properly and are long enough so that no skin is exposed

NOTICE

Risk of equipment damage

Incorrect work practices when handling fluids around the WM-NP3 Mobile Workstation can lead to equipment damage. Observe the following.

- Ensure all surfaces are thoroughly dried before use.
- Use spray type medical agents, such as lubricants, anesthetic or alcohol, away from the WM-NP3 Mobile Workstation.

7.3.2 Cleaning the WM-NP3 Mobile Workstation

Perform the following procedure immediately after use. If cleaning is delayed, organic debris will dry and it can be difficult to effectively clean the surfaces.

1. Take a soft, lint-free wipe that is soaked with a suitable reprocessing agent.
2. Thoroughly wipe the surfaces of the WM-NP3 Mobile Workstation until visibly clean.

3. In case of heavy soiling, use additional clean wipes.
4. Before disinfection, let the surfaces dry or wipe them dry using a lint-free cloth.
5. Visually inspect the product after cleaning. The product must be visually clean.
6. If there are any signs of debris, repeat the cleaning process.

**WARNING****Risk of electric shock**

Using the WM-NP3 Mobile Workstation while it is still wet can result in electric shock. After cleaning, dry the WM-NP3 Mobile Workstation and installed devices thoroughly before use.

7.3.3 Disinfecting the WM-NP3 Mobile Workstation

Low-level disinfection efficacy has been validated by wiping with alkaline disinfectant based on low alcohol (<20%) and quaternary ammonium compounds (Sani-Cloth® Plus manufactured by PDI Professional at 3 minutes contact time).

1. Take a soft, lint-free wipe that is soaked with a suitable disinfectant.
2. Thoroughly wet the surfaces of the WM-NP3 Mobile Workstation.
3. The surfaces must remain visibly wet for the contact time that is specified by the disinfectant's manufacturer.
4. If necessary, use alkaline disinfectant based on low alcohol (<20%) and quaternary ammonium compounds (Sani-Cloth® Plus manufactured by PDI Professional at 3 minutes contact time).
5. Wipe the surfaces dry using a lint-free cloth after the specified contact time.
6. Inspect the products thoroughly before use for traces of wear, as reprocessing can lead to an increased wear of the product.

WARNING

Risk of electric shock

Using the WM-NP3 Mobile Workstation while it is still wet can result in electric shock. After disinfection, dry the WM-NP3 Mobile Workstation and installed devices thoroughly before use.

8 Maintenance

Perform the maintenance described in this chapter to ensure the WM-NP3 Mobile Workstation and accessories remain in a good and safe working condition. Failure to perform maintenance can lead to an unsafe workstation. Observe the following.

- Unless specified otherwise, all repairs should be conducted by Olympus service personnel.
- Maintenance of the separation transformer should only be performed by Olympus service personnel.
- Do not service any parts on the WM-NP3 Mobile Workstation while it is in use with a patient.
- Electrical safety checks should only be performed by suitably qualified personnel.

WARNING

Risk of injury

Should any irregularity be detected during maintenance, do not use the WM-NP3 Mobile Workstation and refer to “Troubleshooting” on page 99 or contact Olympus service center for repair.

Every 6 months

- Check the security and condition of the braked castors. If the braked castors are loose and/or excessively worn, contact Olympus for service.
- Check the power cable and plug for excessive twists, deformation or other irregularities.
- Check the security and condition of the rear panel snaploc couplings. If the snaploc couplings are loose and/or excessively worn, replace with spare parts as detailed in “Customer spare parts” on page 80.

Every 12 months

- Clean each antistatic (conductive) braked castor to remove non-conductive material such as cleaning wax and floor polish, as to not impair the antistatic effectiveness of the braked castors.

- For any faults that occur before the time periods specified for checks, please contact your Olympus service centre.
- Check the screws (and tighten if necessary) on the WM-NP3 Mobile Workstation and any accessories that are fitted. Pay special attention to monitor platforms and the scope poles (horizontal to vertical).

9 Storage and disposal

9.1 Storage

Storing the WM-NP3 Mobile Workstation

1. Disconnect the power cable from the mains electricity.
2. Hang the power cable on the cable hook on the back of the top tray.
3. Move the WM-NP3 Mobile Workstation to the desired storage space.
4. Place the WM-NP3 Mobile Workstation on a level floor.
5. Lock all braked castors.
6. Store the WM-NP3 Mobile Workstation within the environmental conditions that are described in "Technical data" on page 66.

9.2 Disposal

Waste and electric equipment

Refer to your local Olympus distributor for return and/or collection systems available in your country.

Disposal

When disposing of the product, or any of its components, follow all applicable national and local laws, and guidelines.

10 Technical data

10.1 Power cable

WARNING


Risk of injury

To isolate the product from the mains supply, disconnect the mains power plug.

The power cable connects to the separation transformer through a high current IEC-320 style connector, allowing it to be easily replaced if it becomes damaged during use. The power cable is secured to the separation transformer by a pair of clamps, to prevent accidental disconnection during use. To isolate the product from the mains electricity, disconnect the power cable.

10.2 Socket adaptor cable

The socket adaptor cable is designed to provide an interface connection between Olympus endoscopy workstation multiway socket blocks and the power cable supplied with, or integral to, ancillary equipment forming part of the endoscopic system (see illustration).

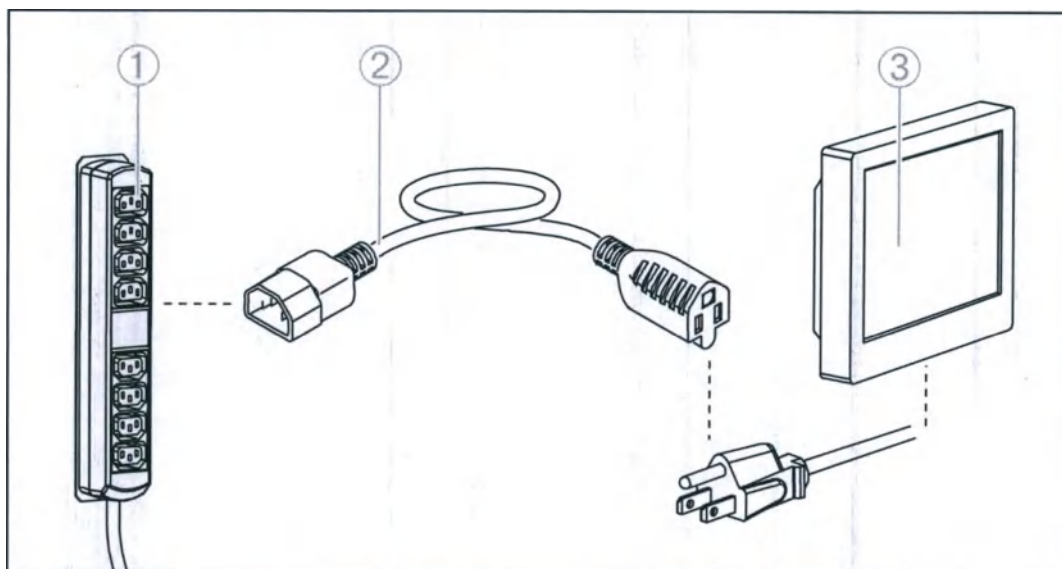- 1) IEC multiway socket block
- 2) Adaptor cable
- 3) LCD monitor

The fully installed system should meet the requirements of IEC 60601-1-1. Refer to "EMC Information" on page 73 for additional information.

The maximum load for the socket adaptor cable is 10 A.

10.3 Protective earth terminal

The protective earth terminal connects to the metal enclosure of the separation transformer and to the earthing system of the output socket.

The line and neutral conductors connect to the primary side of the transformer via two magnetically operated circuit breakers. The circuit breakers have a common actuator and a common inrush limiting circuit.

10.4 Winding

A green LED on the separation transformer illuminates if it is connected to the mains electricity. The separation transformer is wound with a nominal input-to-output ratio of 1:1, with concentric primary and secondary windings. An additional low voltage secondary winding powers the circuit for the power switch of the separation transformer. Electrical separation between windings is maintained by overlapping layers of BoPET (biaxially-oriented polyethylene terephthalate).

10.5 Primary winding

The primary winding draws an inrush current surge when it is connected to the mains electricity, depending on the magnetic state of its core at the instant when it is powered. The inrush limiting circuit inserts a high power thermistor in series with the primary winding to limit the initial surge, dissipating it as heat. After approximately 0.3 seconds, a relay closes (producing an audible 'click') to by-pass the thermistor for the following reasons.

- To avoid steady-state losses that would adversely affect the load regulation.
- To allow the thermistor to cool for the next switch-on.

10.6 Secondary winding


The secondary winding supplies power to the separated output-sockets. An additional output socket connects to the power switch located on the front of the workstation's top tray, operating a relay to control the power to the output-sockets.

10.7 Potential equalization terminals

A potential equalization terminal is provided next to the input cable for use where local regulations require the enclosures of each part of the system be connected to a common earth point, independently of the power cables. Additional terminals can be installed by attaching a six-way equipotential terminal strip, available as an accessory, to each side of the separation transformer.

10.8 Block diagram

The block diagram shown shows the different sub-units within the separation transformer.

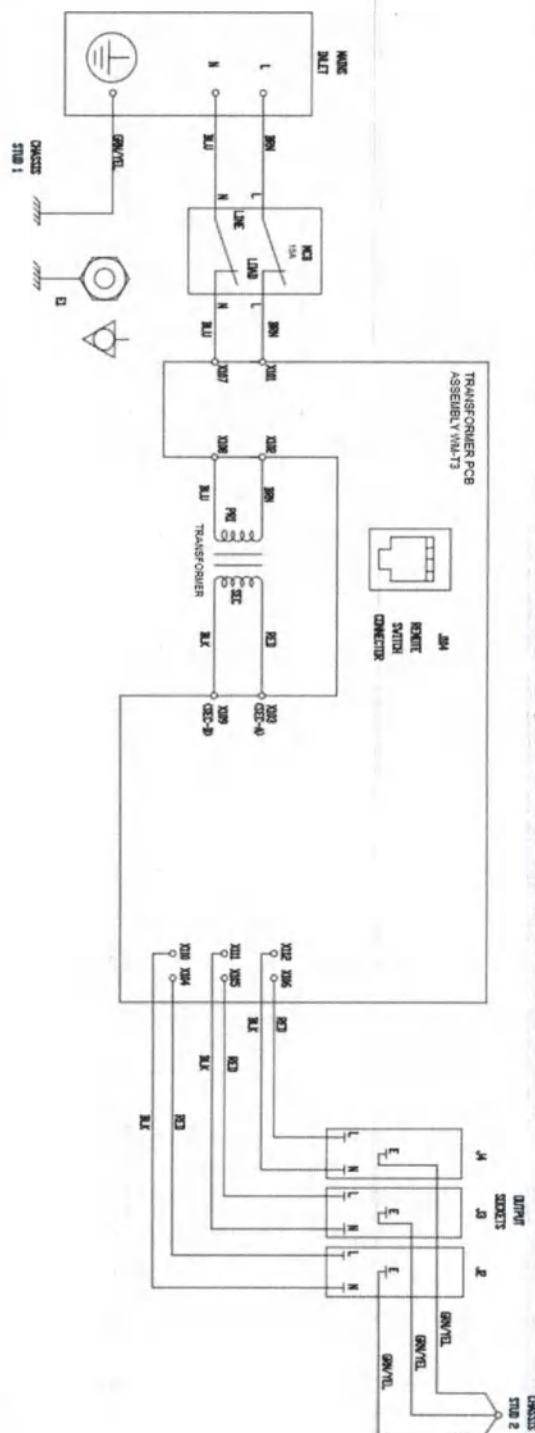10.9 Safety standards for medical electric systems

CAUTION

Risk of injury

Failure to observe the information in this section can result in the separation of the entire system being compromised, defeating the protection provided by the transformer in a single fault condition

The following are standards applicable to medical electrical equipment.

- EN 60601-1:2006+A12:2014
- IEC 60601-1:2005
- UL 60601-1:2003
- CAN/CSA C22.2 601.1 M90

The following are standards applicable to medical electrical systems.

- Clause 16 of EN 60601-1:2006+A12:2014
- IEC 60601-1-1:2005

Earth leakage current of medical electrical systems, ancillary equipment and medical electrical equipment

Regulatory standards require that medical electrical systems meet the earth leakage current requirements of EN IEC 60601-1/UL 60601-1, which allows a maximum earth leakage current of 0.5 mA in normal condition (0.3 mA in North America). Combined earth leakage current can result in excessive enclosure leakage current under single fault conditions, such as the interruption of a common protective earth path.

A separation transformer is required under one of the following conditions.

- Use of ancillary equipment that does not meet the earth leakage current requirements within the 'patient environment'. The separation transformer supplies ancillary equipment and medical electrical equipment, as long as the maximum separation transformer rating is not exceeded.
- The total earth leakage current of the medical electrical system is above 0.5 mA (0.3 mA in North America). The separation transformer must be used to supply the medical electrical equipment sufficiently, therefore reducing the total earth leakage current to below the allowed limit, as long as the maximum separation transformer rating is not exceeded.

NOTE

The total impedance of the protective earth path for a system can be up to 0.4 W or higher if the conditions of clause 8.6.4b) of EN 60601-1 / 18 g) of IEC/UL 60601-1 are satisfied.

Power all medical electrical equipment from the separation transformer, unless appropriate separation is employed as described in Annex I of EN IEC 60601-1:2006 / Annex EEE of IEC 60601-1-1:2005. This particularly applies where ancillary equipment is remotely situated, or where signal input and signal output parts are coupled between the medical electrical equipment and a facility's data network.

Use of high Intensity light sources

High intensity light sources use arc lamps, which require a high inrush current. Do not supply high intensity light sources from the the separation transformer if the maximum separation transformer rating is exceeded.

Use of multiple separation transformer

Some system installations, for instance those with 2 mobile workstations, require 2 separation transformers to be fitted. Do not connect extension cables and additional multiple socket outlets to the system. Only use the output-sockets of the separation transformer to supply medical electrical equipment that forms part of the endoscopy imaging system.

10.10 EMC Information

10.10.1 EMC (electromagnetic compatibility)

The WM-NP3 Mobile Workstation is EMC-tested in conformity with the requirements of IEC 60601-1-2 Edition 4 and can be used in the vicinity of other EMC-tested devices that fulfill the requirements of the relevant IEC 60601-1-2 Edition 4 standard. The WM-NP3 Mobile Workstation is a medical device that requires special safety precautions and must be installed and placed in operation according to the following EMC information on portable and mobile RF communication devices as detailed in this chapter.

10.10.2 Essential performance with regards to EMC

Olympus has determined that the WM-NP3 Mobile Workstation does not have any functions, which if they failed to work properly, would impact the safety of the system. All necessary instruction for maintaining essential performance with regard to electromagnetic disturbances for the expected service life of the unit is listed in this chapter.

NOTE

Output power of the transformer may be reduced for a period of 120 ms following short dips in the mains electricity.

10.10.3 Environment with regards to EMC

The WM-NP3 Mobile Workstation is suitable for the professional health care facility environment as defined in IEC 60601-1-2 Edition 4. In addition to this, the WM-NP3 Mobile Workstation can operate near Olympus endoscopic HF surgical equipment outside a RF shielded room. Therefore, for this purpose only, the WM-NP3 Mobile Workstation is suitable for the special environment as defined IEC 60601-1-2 Edition 4. For maintaining basic safety and essential performance with regard to electromagnetic disturbances use only Olympus specified parts.

WARNING
Risk of injury

Use of equipment adjacent to or stacked with non-Olympus equipment should be avoided because it can result in improper operation. If such use is necessary, observe the equipment to verify that it is operating normally.

WARNING
Risk of injury

Portable RF communications equipment (including peripherals such as antenna cables and external antennas) should be used no closer than 30 cm (12 inches) to any part of the WM-NP3 Mobile Workstation including cables specified by the manufacturer. Otherwise, degradation of the performance of this equipment could result.

10.10.4 Electromagnetic emissions

The WM-NP3 Mobile Workstation is intended for use in the electromagnetic environment specified below. Ensure that the WM-NP3 Mobile Workstation is used in such environment.

Emission tests compliance	Level
RF emissions CISPR 11	Group 1 Class B
Harmonic emissions IEC 61000-3-2	Class A
Voltage fluctuations/flicker emissions IEC 60000-3-3	Complies

10.10.5 Electromagnetic immunity

The WM-NP3 Mobile Workstation is intended for use in the electromagnetic environment specified below. Ensure that the WM-NP3 Mobile Workstation is used in such environment.

Immunity	IEC 60601-1-2 test level
Electrostatic Discharge (ESD) IEC 61000-4-2	± 8 kV contact, indirect ± 15 kV air
Electrical fast transient/burst IEC 61000-4-4	± 2 kV for power supply lines ± 1 kV for input/output lines
Surge IEC 61000-4-5	± 1 kV line(s) to line(s) ± 2 kV line(s) to earth
Voltage dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11	1) 0% U_n 1 cycle 2) 0% U_n ½ cycle, Phase 0,45,90,135,180,225,270,315 3) 70% U_n (30% dip on U_T for 25 cycles) 4) 0% U_T for 5 sec.
Conducted RF IEC 61000-4-6	6 Vrms ISM bands 3 Vrms 150 kHz to 80 MHz
Radiated RF IEC 61000-4-3	3 V/m 80 MHz to 2.7 GHz IEC 60601-1-2:Ed 4 Table 9 bands/levels/modulations.

10.10.6 Accessories, transducers and cables

WARNING

Risk of injury

Use of accessories, transducers and cables other than those specified in this document can result in increased electromagnetic emissions or decreased electromagnetic immunity of this equipment and result in improper operation.

10.11 Specifications

10.11.1 WM-NP3 Mobile Workstation

The following table describes specifications for the WM-NP3 Mobile Workstation.

Product name	WM-NP3 Mobile Workstation	
Dimensions	Height	2140 mm (scope pole in highest position)
		1400 mm (workstation only)
		1150 mm (to top tray)
		2073 mm (with 26" LCD monitor in highest position)
	Depth	675 mm
	Width	665 mm
Weight	84 kg unladen, including fitted separation transformer	
Maximum weight capacity (evenly distributed)	311 kg	
Load capacity (evenly distributed)	Top tray	20 kg
	Intermediate shelves	31 kg
	Base panel	35 kg
	NOTE: Equipment stored on the WM-NP3 Mobile workstation should not exceed the size of the shelf.	
Castors	4 x 125 mm twin wheel, 2 x antistatic (conductive) brakes and 2 x non-antistatic (non-conductive) brakes.	
Ambient temperature	Storage	-40 °C to +70 °C (-40 °F to +158 °F)
	Operation	+10 to +40 °C (+50 to +104 °F)
Relative humidity	95% maximum at +40 °C (+104 °F) non-condensing	
Atmospheric pressure	Storage / Operation	70 to 106 kPa (10.15 to 15.37lb/in ²)

10.11.2 Separation transformer

The following table describes specifications for the separation transformer.

Product name	Olympus separation transformer				
Identification	Market	Japan	USA		UK, Europe, Rest of the world
	Model no.	MAJ-2155 (PSE compliant)	MAJ-2156	MAJ-2156-15A	MAJ-2157
	Part no.	K10028136	K10028137	K10037151	K10028138
Power input requirements	Input voltage	100 V	110-120 V	110-120 V	220-240 V
	Frequency	50/60 Hz	60 Hz	60 Hz	50/60 Hz
	Power input (max)	1400 W	1650 VA	1300 VA	1900 VA
	Cable part no.	K10001001	K10034395	K10001002	K10001003 (UK) K10001004 (Europe) K10001005 (RoW)
Inrush current	At maximum rated input voltage, cold start at 25 °C ambient temperature	25 A	25 A	25 A	25 A
Power output	Maximum load	1300 VA	1550 VA	1250 VA	1800 VA
	Maximum load current	14 A (Maximum: 10 A per outlet-socket)	16 A	12 A	8 A
Overload protection	Circuit breakers (dual circuit, single actuator) breakers (dual circuit, single actuator)	2 x 15 A	2 x 16 A	2 x 13 A	2 x 9 A
Potential equalization	Terminals fitted	1 at input			

Classification	Standards compliance	The separation transformer conforms to the relevant requirements of: EN 60601-1, IEC 60601-1, ANSI/AAMI ES 60601-1, UL 60601-1, CAN/CSA C22.2 No. 60601-1, CAN/CSA C22.2-601.1
	Type of protection against electric shock	Class I equipment
	Degree of protection against explosion	None, the separation transformer should not be used within the zone of risk of flammable gases
	Degree of protection against ingress of water	IPX0 as defined in EN IEC 60529:1992, ordinary equipment
	Mode of operation	Continuous
Environmental conditions	Ambient temperature	Storage: -40 °C to +70 °C (-40 °F to +158 °F)
		Operation: +10 to +40 °C (+50 to +104 °F)
	Atmospheric pressure	Storage/Operation: 70 to 106 kPa (10.15 to 15.37 lbf/in ²)
	Relative humidity	95% maximum at +40 °C (+104 °F) non-condensing
Year of manufacture	 The year of manufacture is shown on the packaging label, combined with the 'MANUFACTURER' symbol.	
End of life	 The product must not be disposed of as unsorted municipal waste, but should be collected separately. Refer to your local Olympus distributor for return and/or collection systems available in your country.	

11 Customer spare parts


The following parts have been made available as replacement parts. Please contact Olympus service center to order any required parts. In order to maintain standards compliance, it is important that replacement components should have the same specification as those fitted originally. This can best be guaranteed by using only those spare parts and components listed in this instruction for use.

Only those spare parts listed below should be fitted by the customer.


11.1 WM-NP3 Mobile Workstation spares

Spare part	Part number
Base bumper	K10029937
Water bottle bracket	K10036097
Fitting kit	K10036057
Instructions for use	K10035969
Power cable (UK)	K10028273
Power cable (Flying Leads – no Plug)	K10028275
Power cable (USA) - 20 A	K10034852
Power cable (USA) - 15 A	K10028272
Power cable (Japan)	K10035719
Power cable (Schuko)	K10028274
IEC cable (0.8 m)	K10028276
IEC cable (1.5 m)	K10029375
IEC cable (2.5 m)	K10000449
Rear panel	K10036095
Snaploc coupling 8 mm	K10036096
Cable management	K10036094
Socket adaptor cable (USA)	K10036766
Socket adaptor cable (Japan)	K10036760


Replacing the base bumper

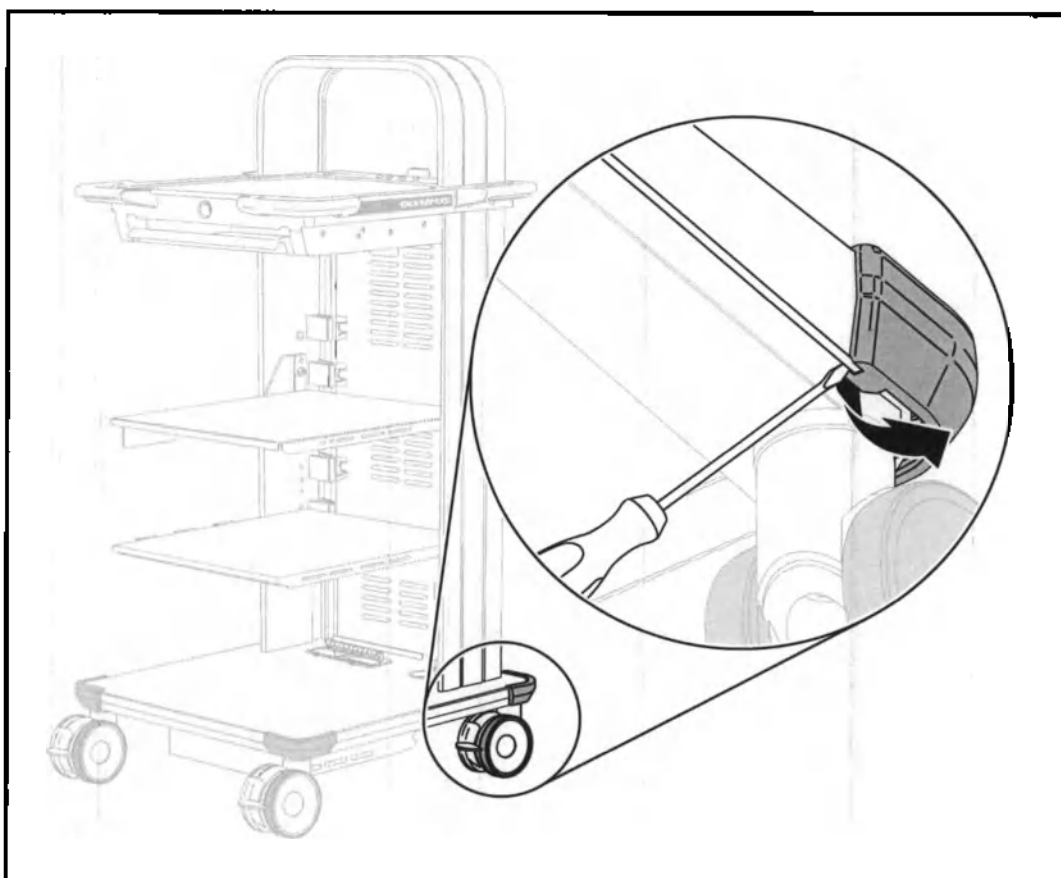1. From the underside of the base panel, loosen the 2 hooks on the base bumper from the slots in the base panel with a flat-blade screwdriver.

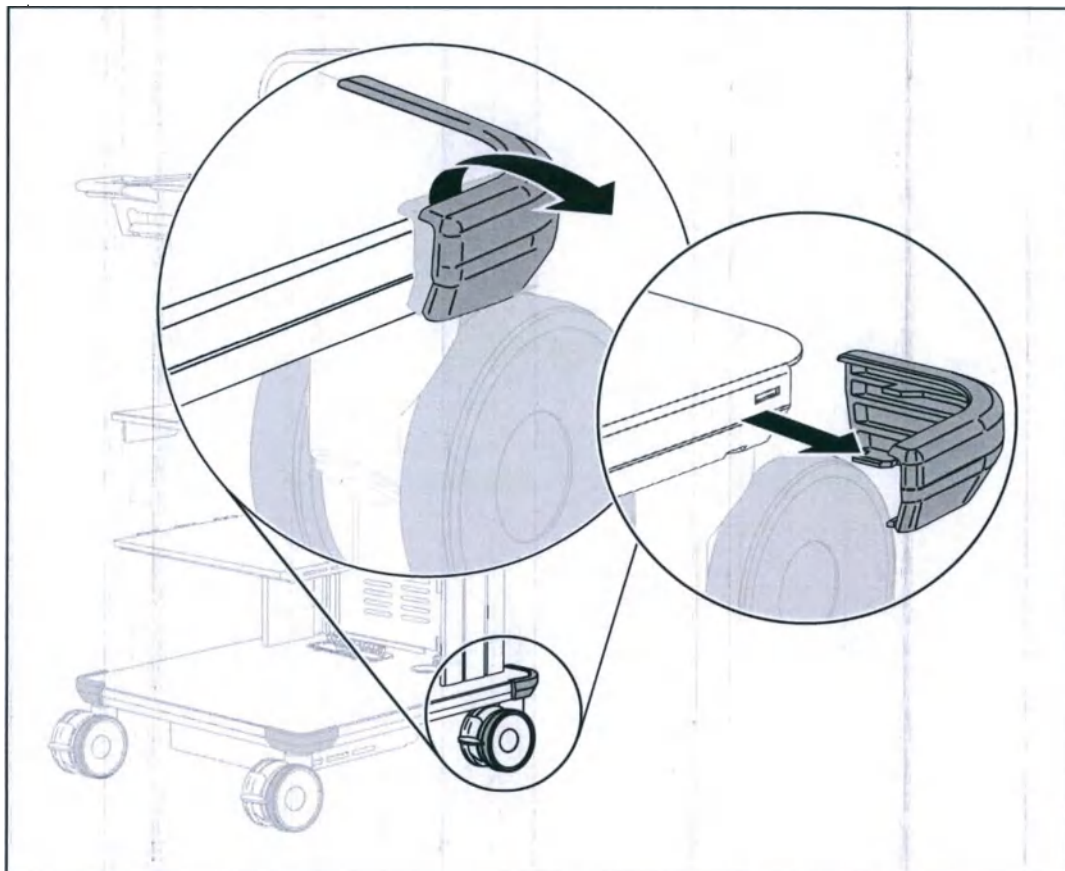
2. Pull the base bumper away.

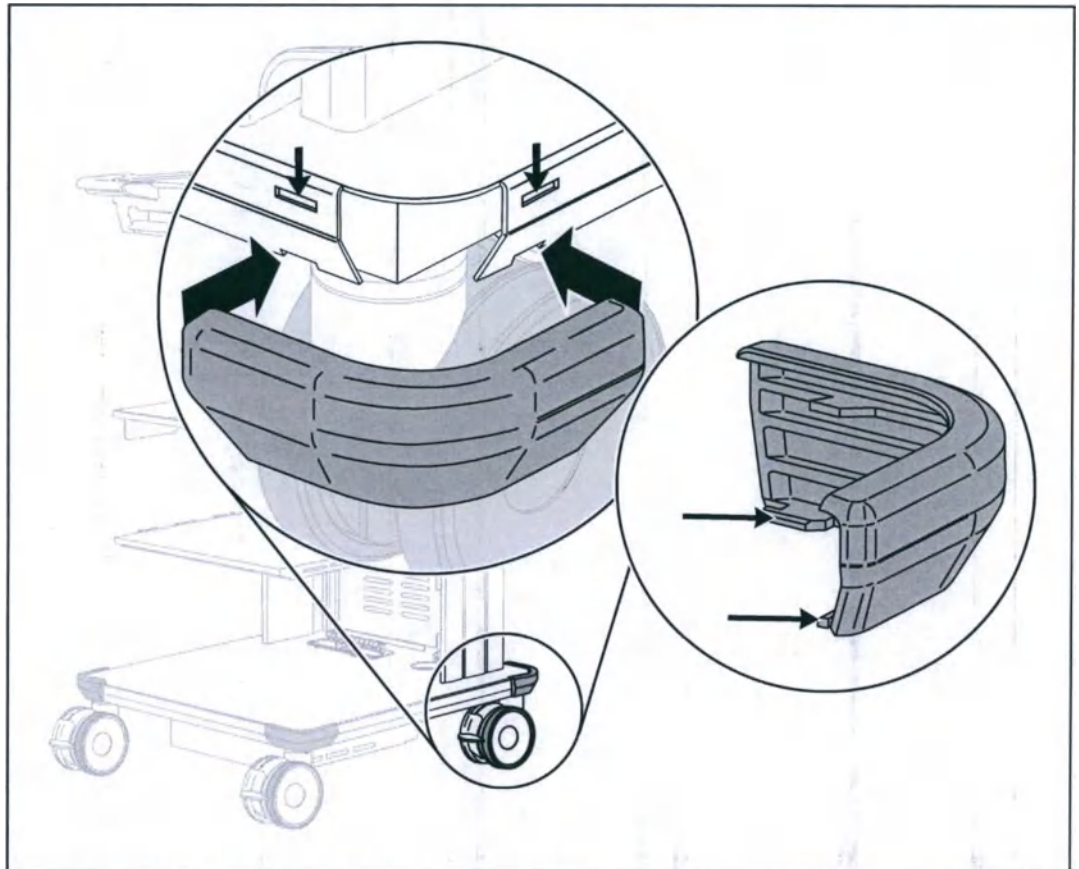
3. Fit the replacement base bumper by aligning the 2 internal hooks with the slots on the adjacent side faces of the base panel.
4. Push the base bumper from the corner until both lower hooks have engaged with the slots on the underside of the base panel (audible click).
5. Confirm the base bumper is fully seated before use.


Replacing the power cable


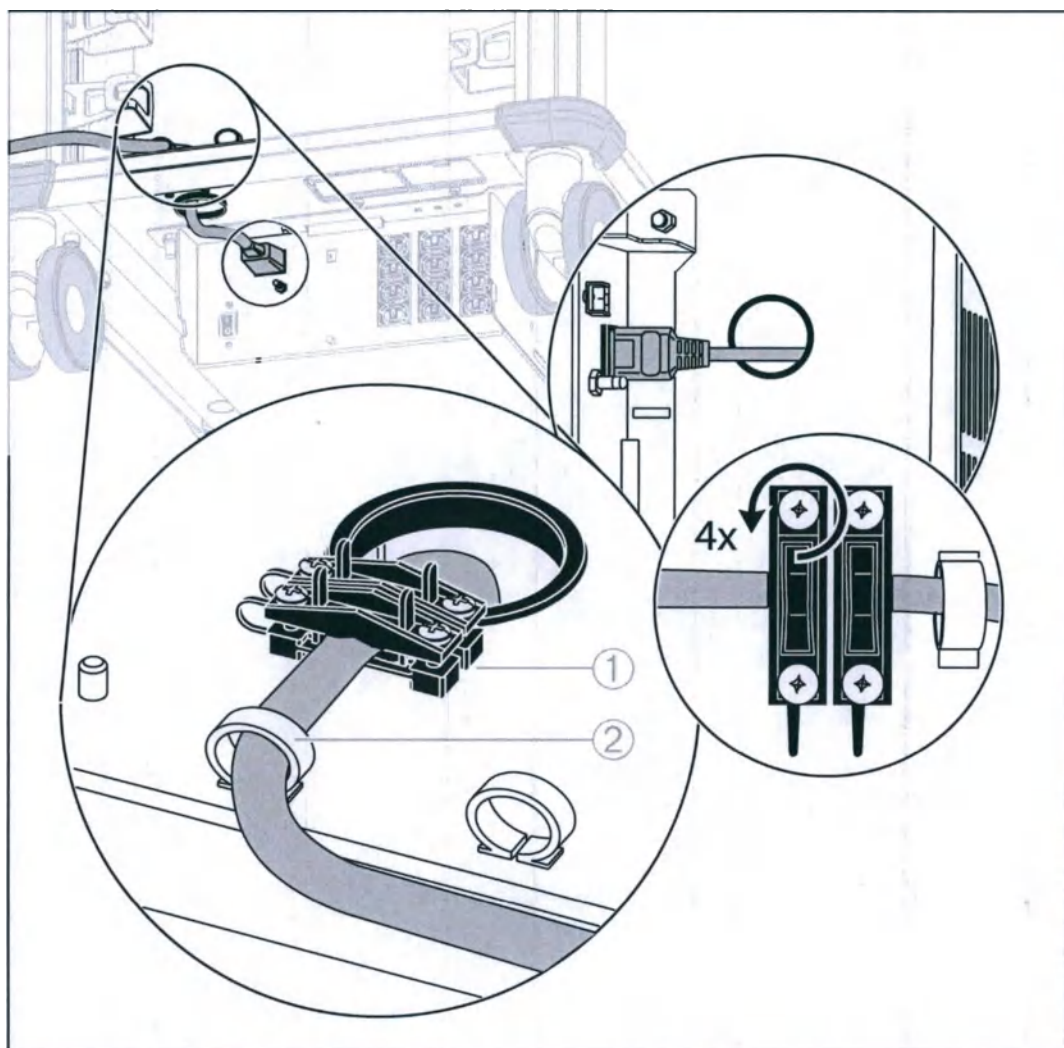
- 1) Black clamp
- 2) Power cable clip
- 1. Loosen the 4 screws from the 2 black clamps retaining the power cable at the rear of the base panel.

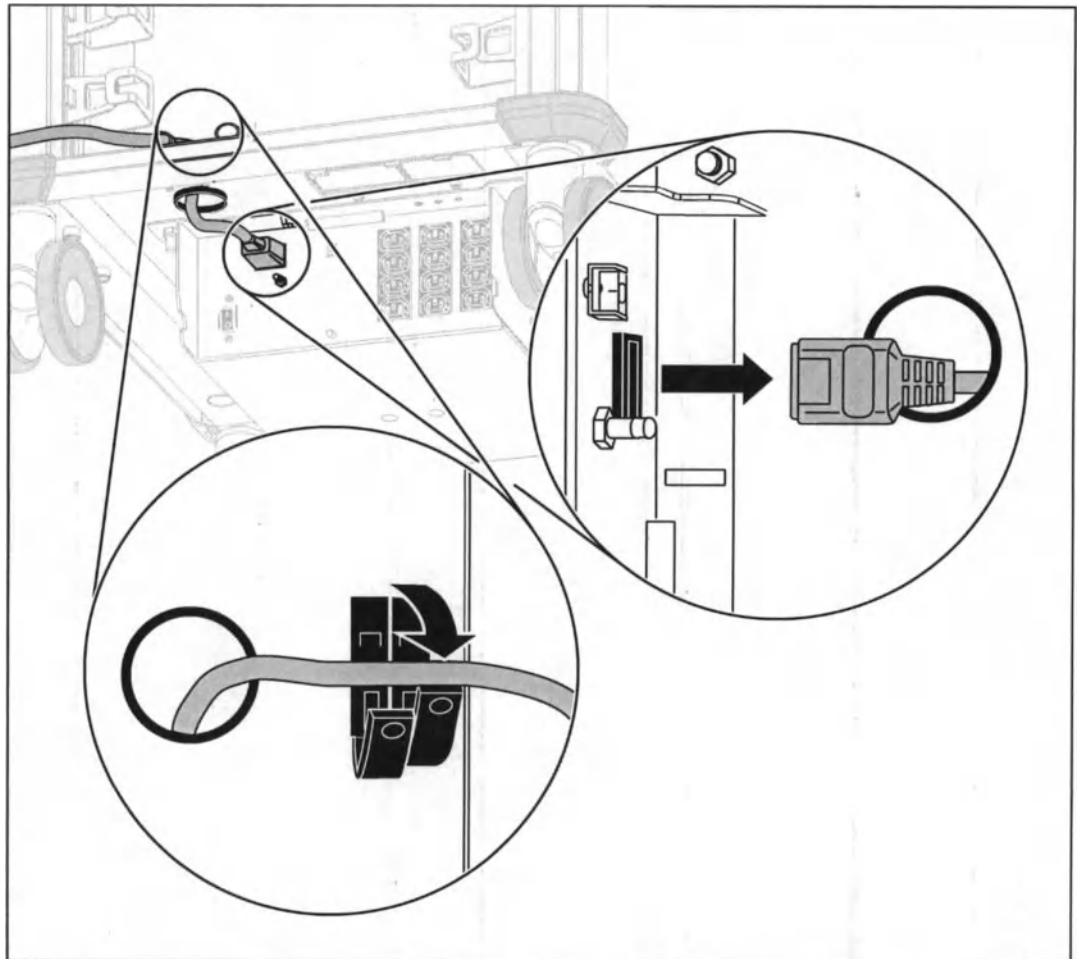
2. Open the clamps and disconnect the power cable from the separation transformer.

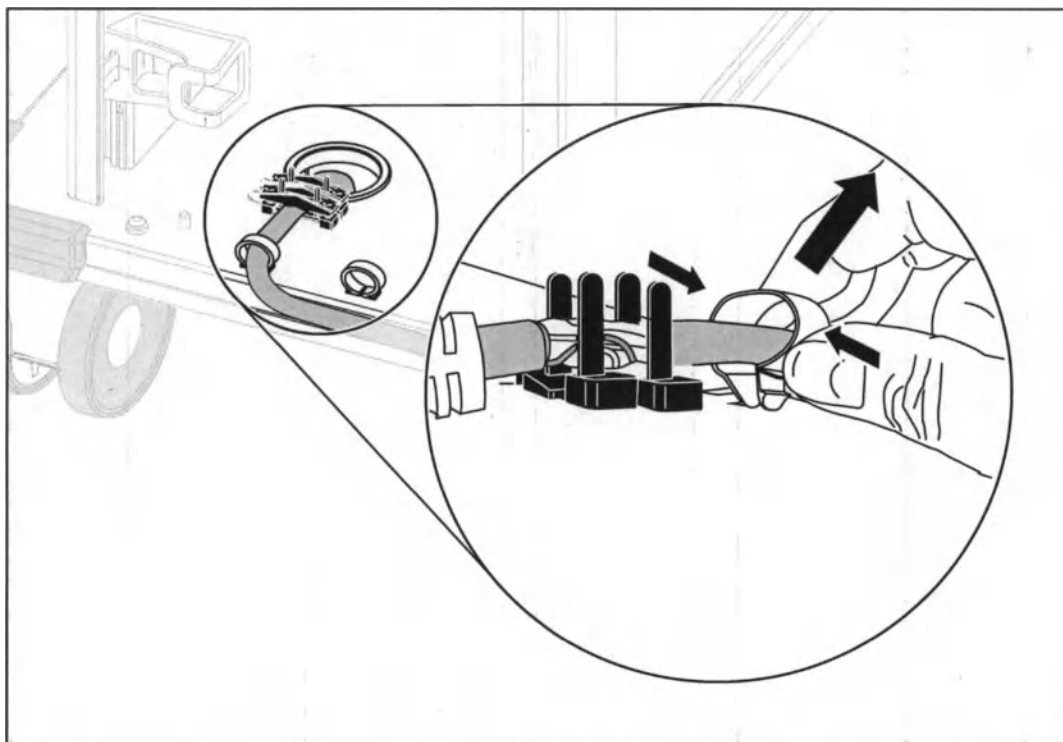
3. From the upper side of the base panel, locate the position of the power cable clip and squeeze the prongs together, releasing the cable clip from the base panel.

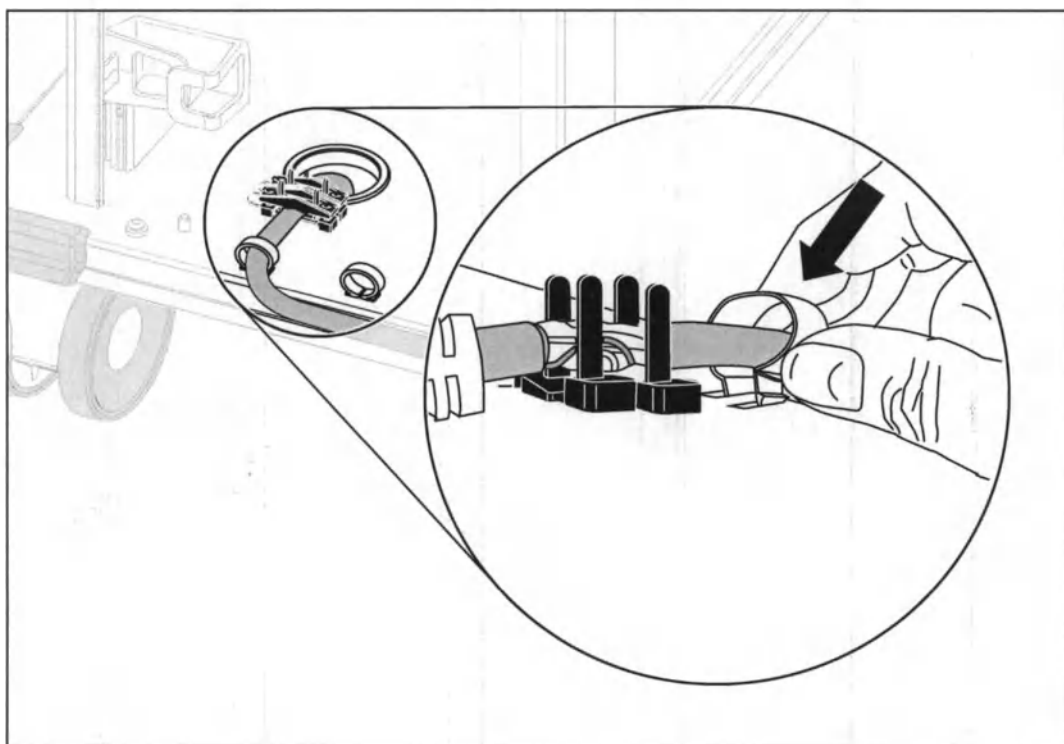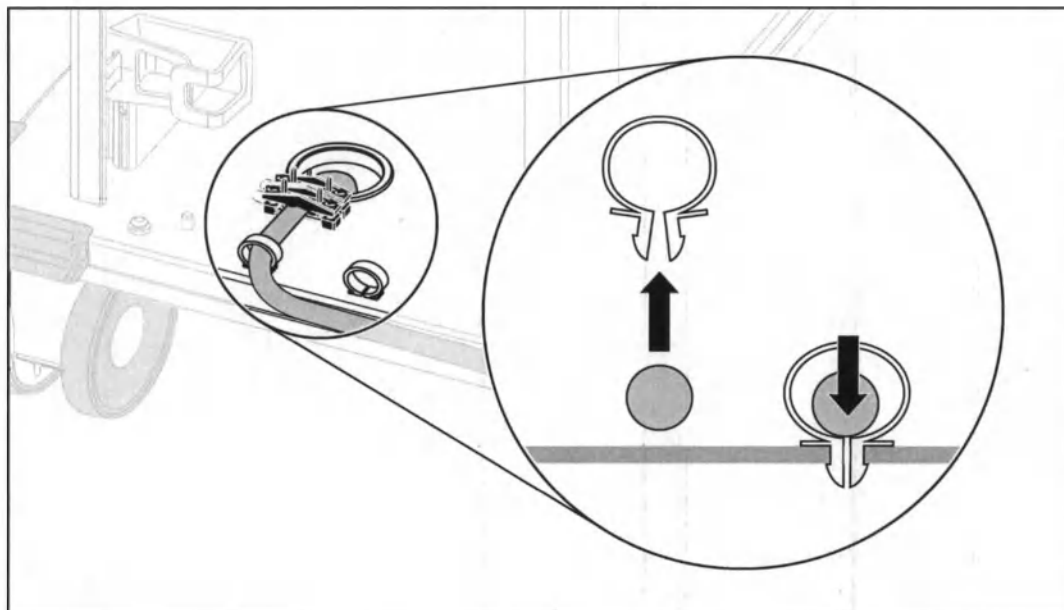4. Place the new power cable through the cable clip and push the cable clip prongs into the base panel hole, encapsulating the power cable.

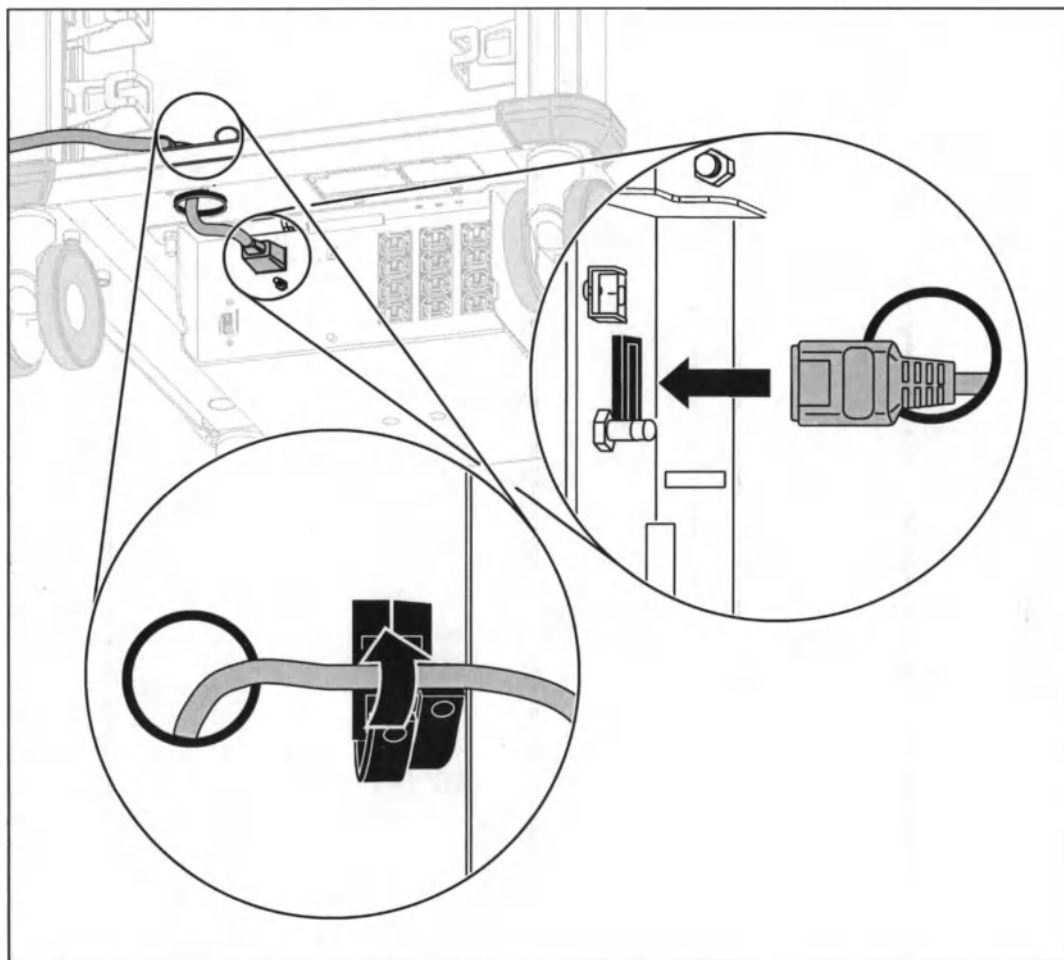5. Push the power cable through the large round aperture and connect to the connector on the separation transformer, ensuring that the connectors are fully seated.

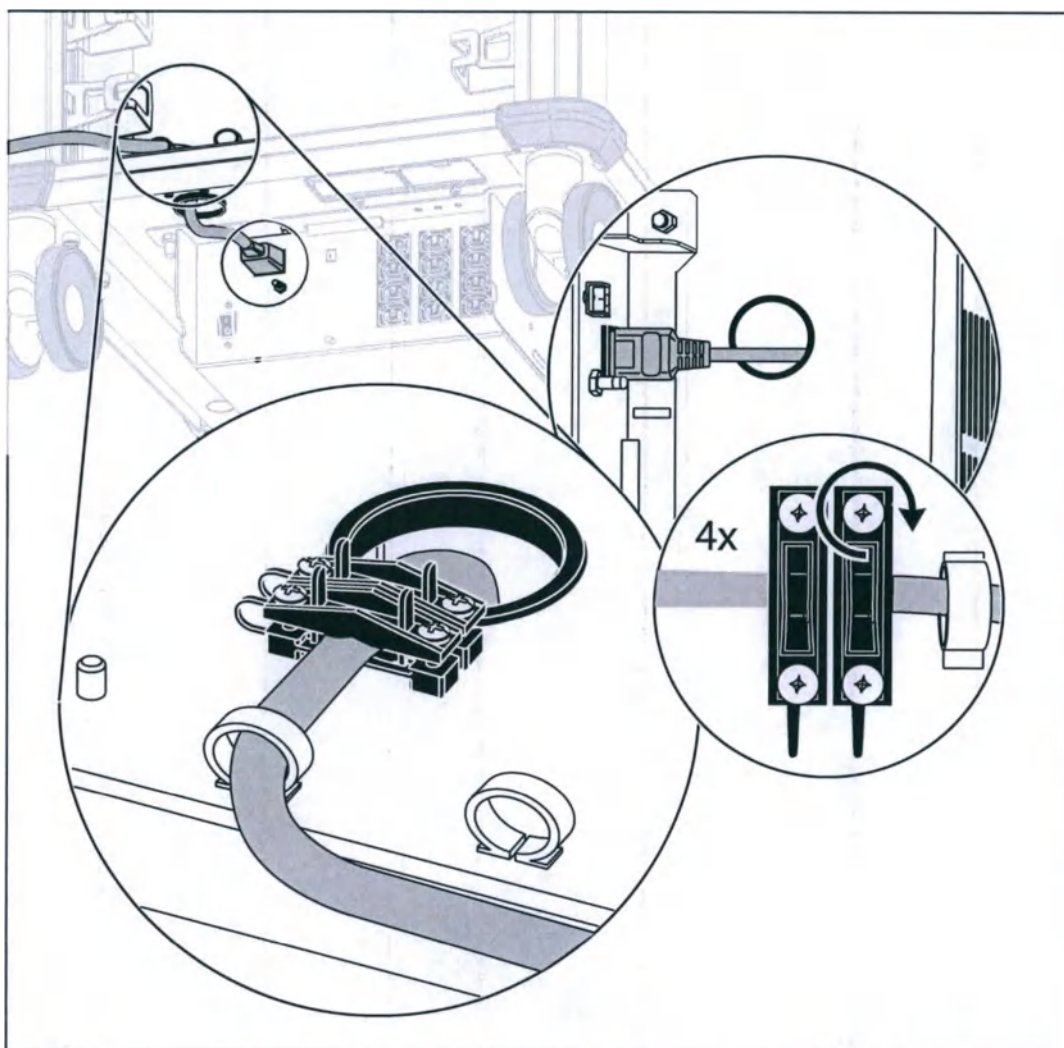
6. Lay the power cable within the 2 black clamps and secure the clamps with the 4 screws, ensuring that the power cable is fully restrained and that the connection to the separation transformer has been maintained.

Replacing the IEC cable

1. Carefully cut any cable ties used to position or restrain the IEC cable to be replaced.
2. Disconnect the IEC cable from both the medical electrical equipment and the separation transformer.
3. Carefully remove the IEC cable from the workstation cable management so as not to disturb existing cabling.


4. Connect the replacement IEC cable to both the medical electrical equipment and the separation transformer, securing the cable to the strain relief moulding on the base panel with a cable tie.

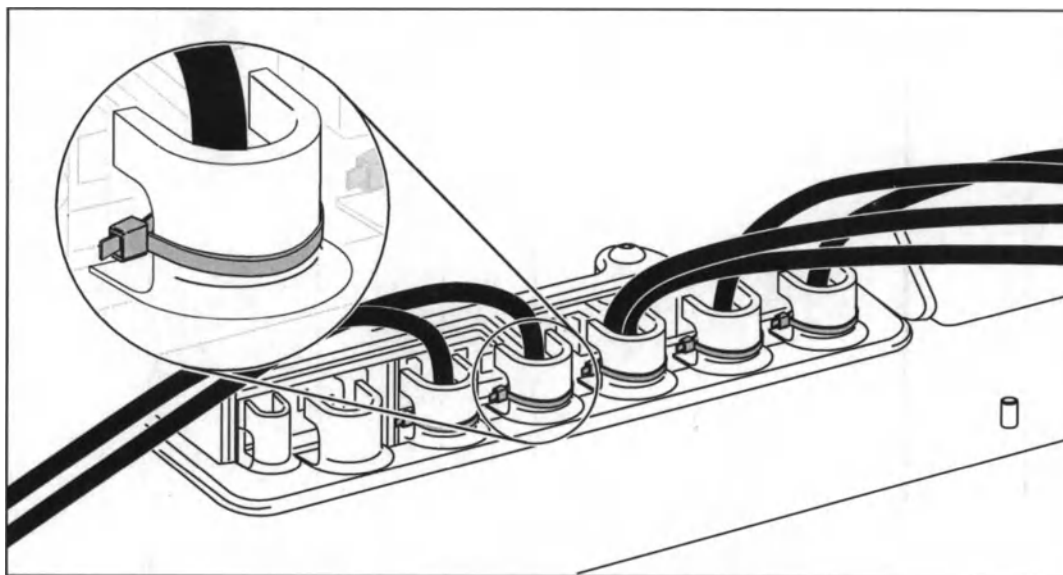
5. Dress the cable within the trolley cable management and re-secure.

Replacing the cable management

1. Remove the rear panel, refer to "Removing and attaching the rear panel" on page 30.

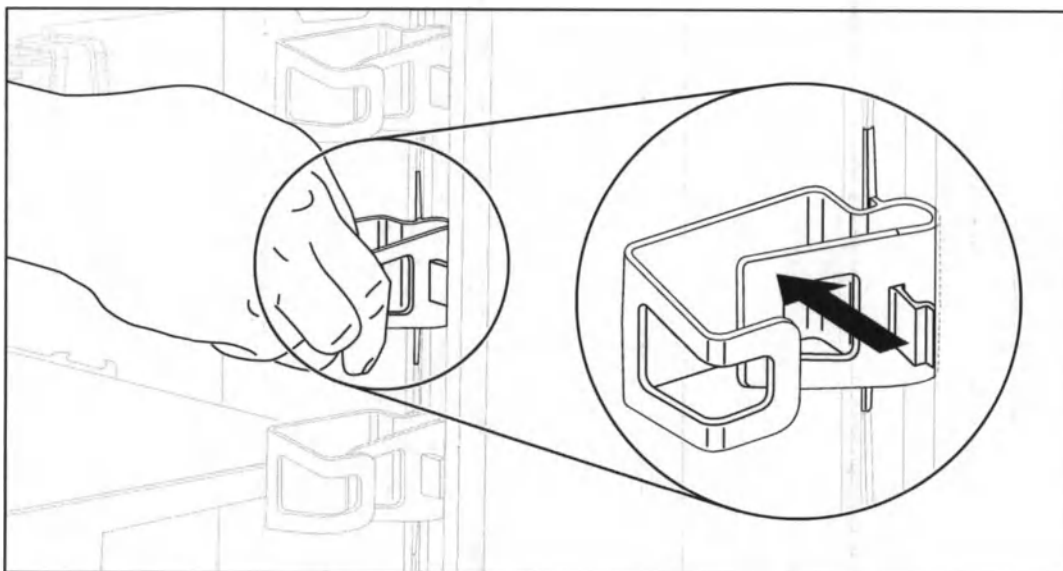2. Push one unit of the cable management together.

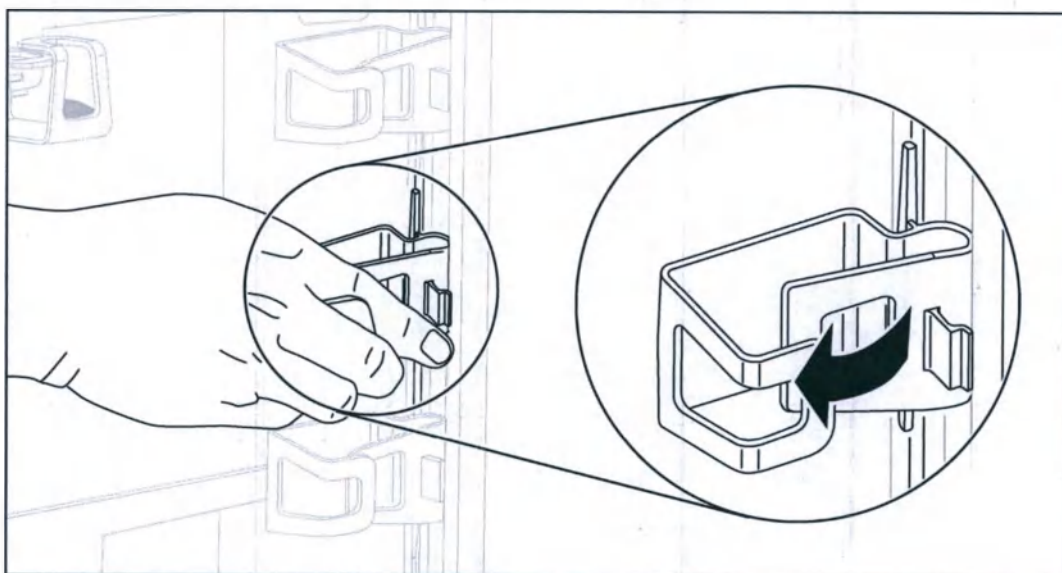
3. Tilt the cable management unit inwards to remove it from the workstation hoop.
4. Push the new unit of the cable management together and place it in the slot.
5. Let go of the new unit of the cable management.

11.2 Scope pole kit spares


Spare part	Part number
Endoscope support	K10029949
Endoscope hanger assembly	K10029954
Fitting kit	K10036060
Instructions for use	K10035974

Replacing the endoscope support

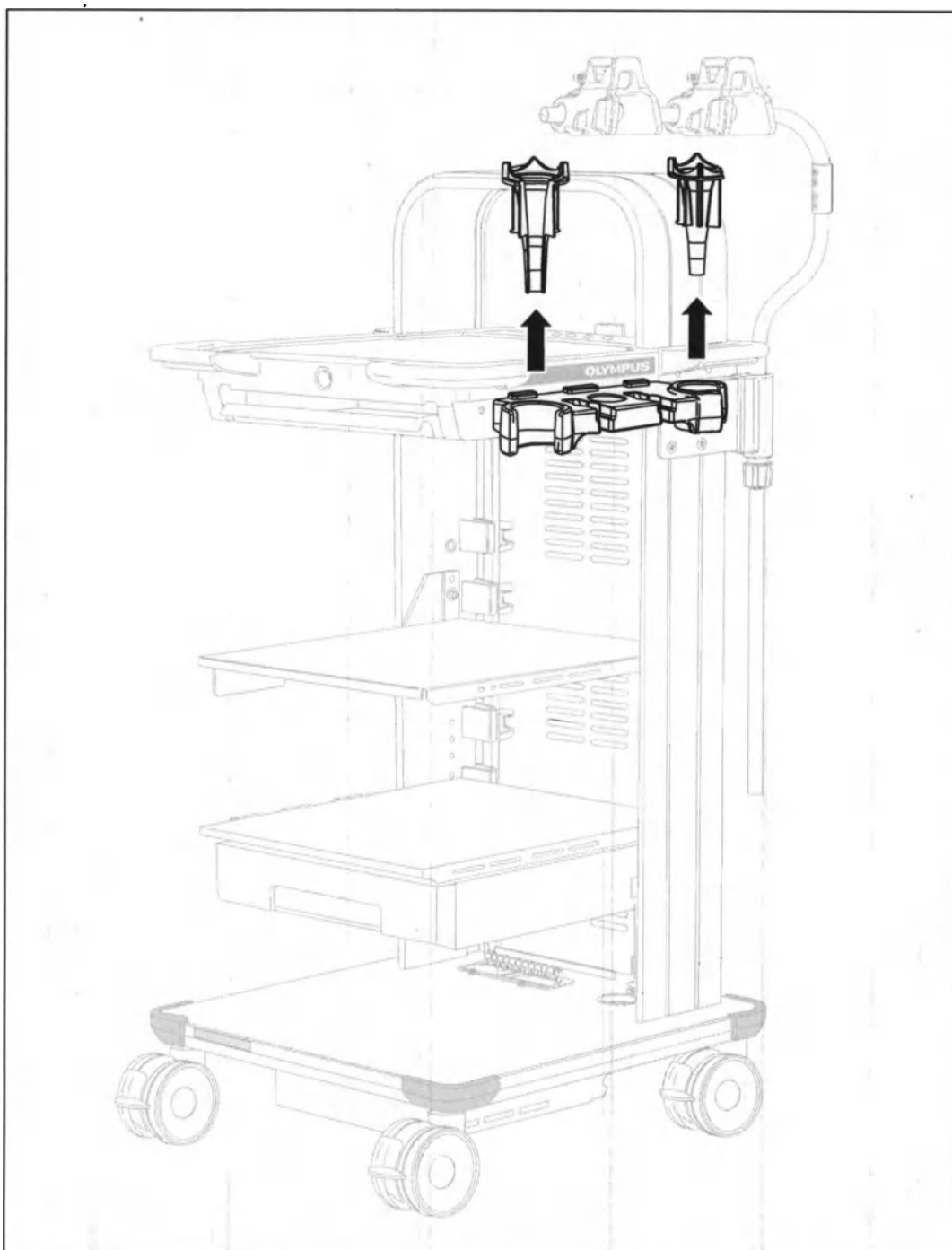
1. Squeeze the damaged endoscope support together and remove them from the endoscope connector holder.
2. Squeeze the replacement endoscope support together and insert them into the endoscope connector holder.
3. Ensure the endoscope supports are fully seated prior to use.

Replacing the endoscope hanger

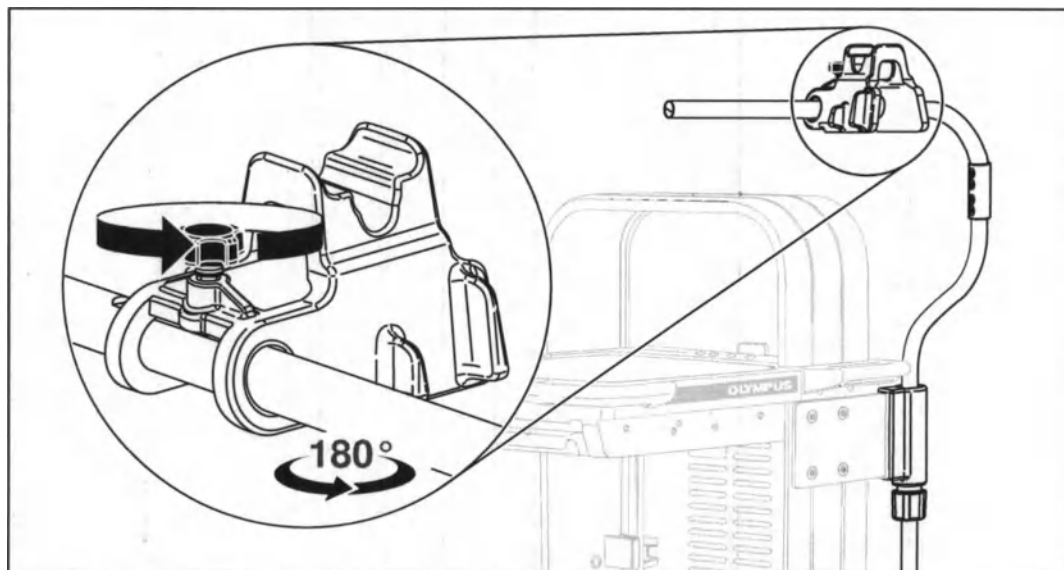- 1. Release the handwheel on the top of the endoscope hanger.**

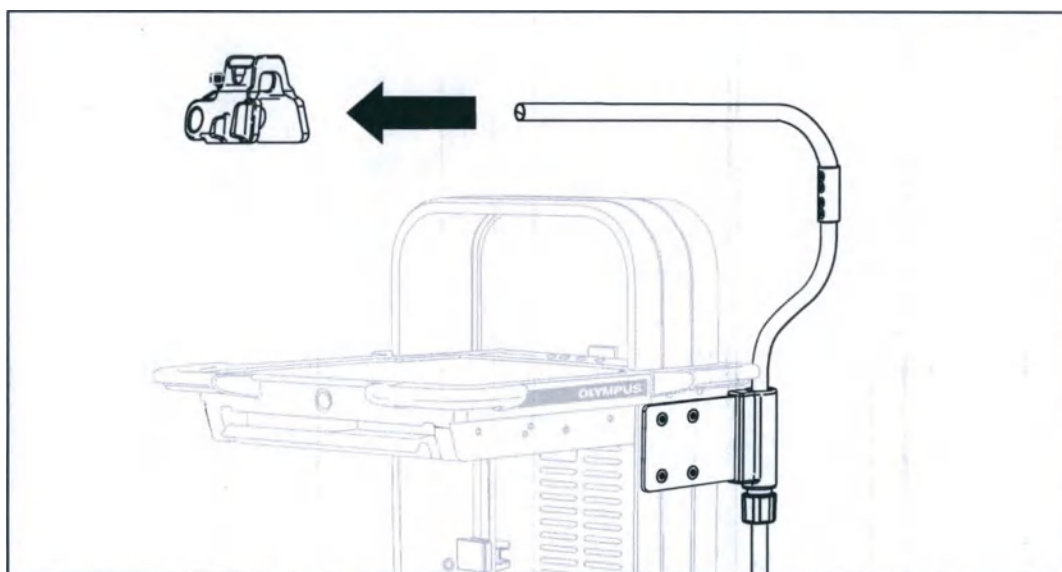
2. Slide the endoscope hanger off the endoscope pole.
- 3. Slide the replacement endoscope hanger onto the endoscope pole.**
- 4. Tighten the handwheel of the endoscope hanger to secure.**


11.3 Sliding keyboard tray spares

Spare part	Part number
Snaploc coupling	K10029941
Fitting kit	K10029942
Instructions for use	K10029943

Replacing the snaploc coupling

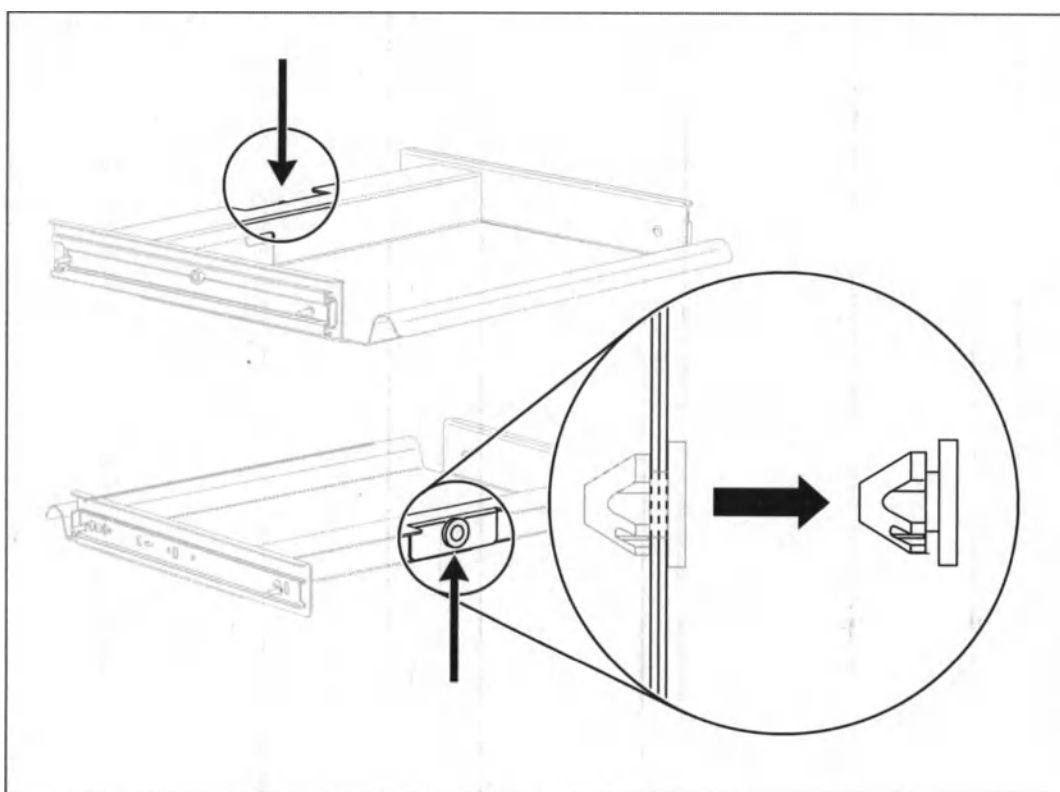1. Loosen the snaploc couplings from the keyboard tray aperture with your fingers, and discard.

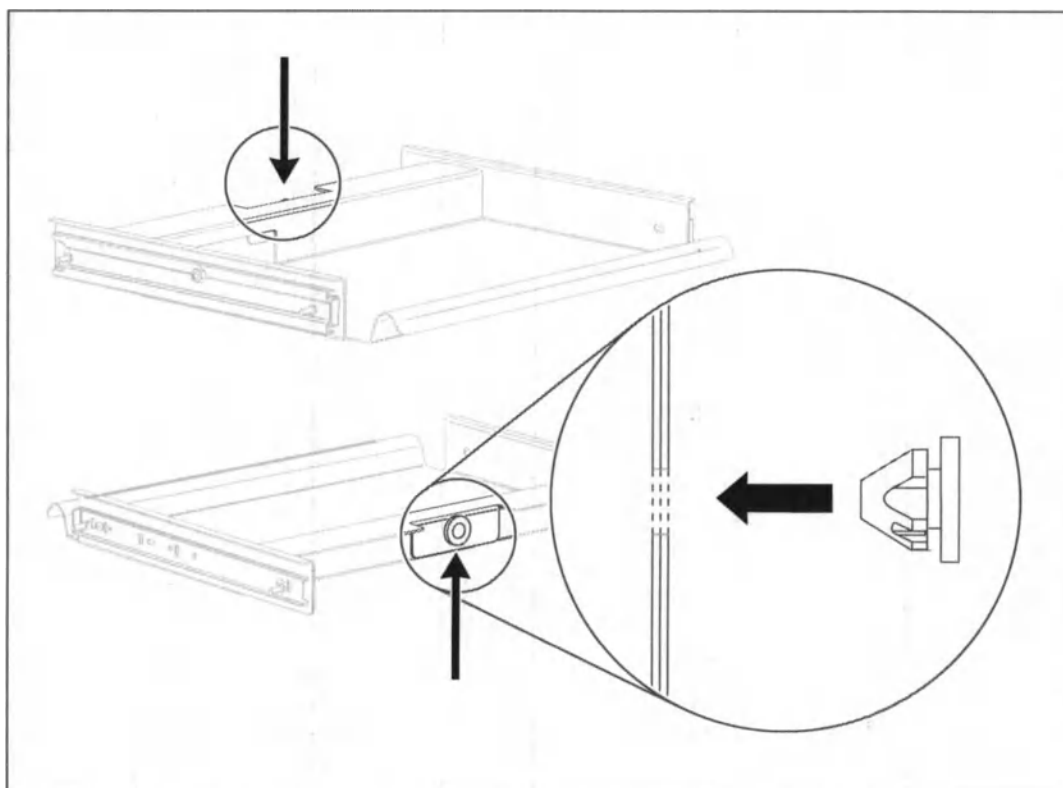
2. Insert the new snaploc coupling into the hole and press firmly until retained.


11.4 Rear panel spares

Spare part	Part number
Snaploc coupling 8 mm (pair)	K10036096

Replacing the snaploc coupling

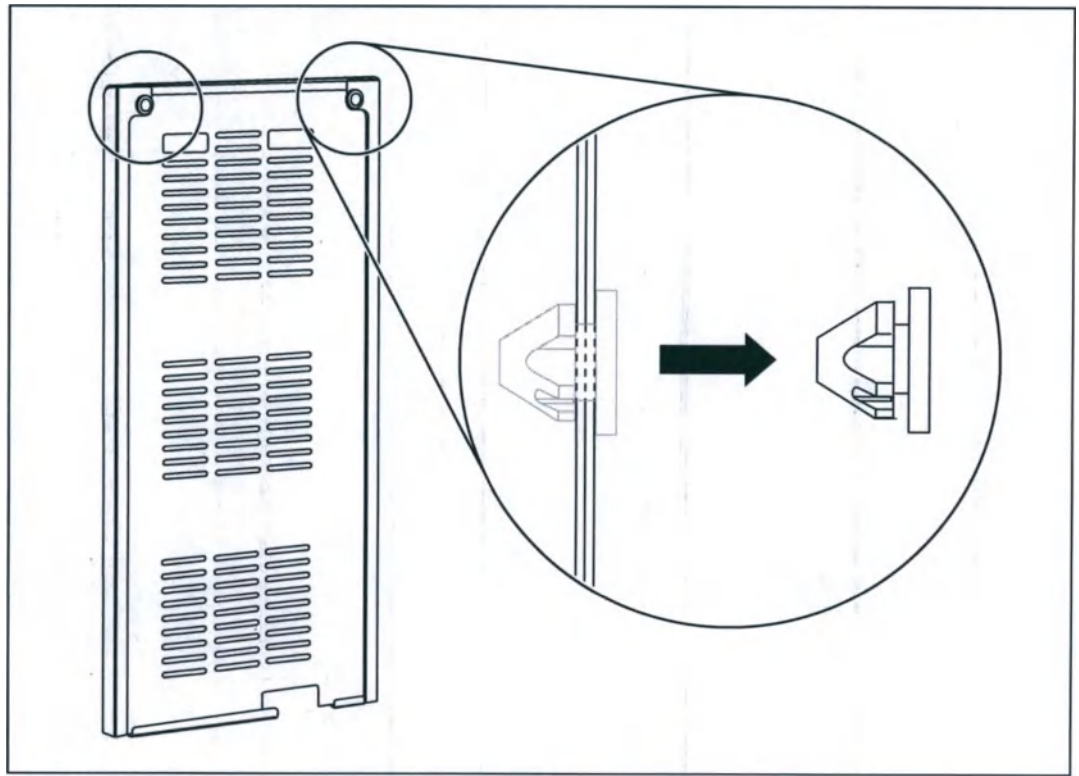1. Loosen the snaploc couplings from the rear panel with your fingers, and discard.

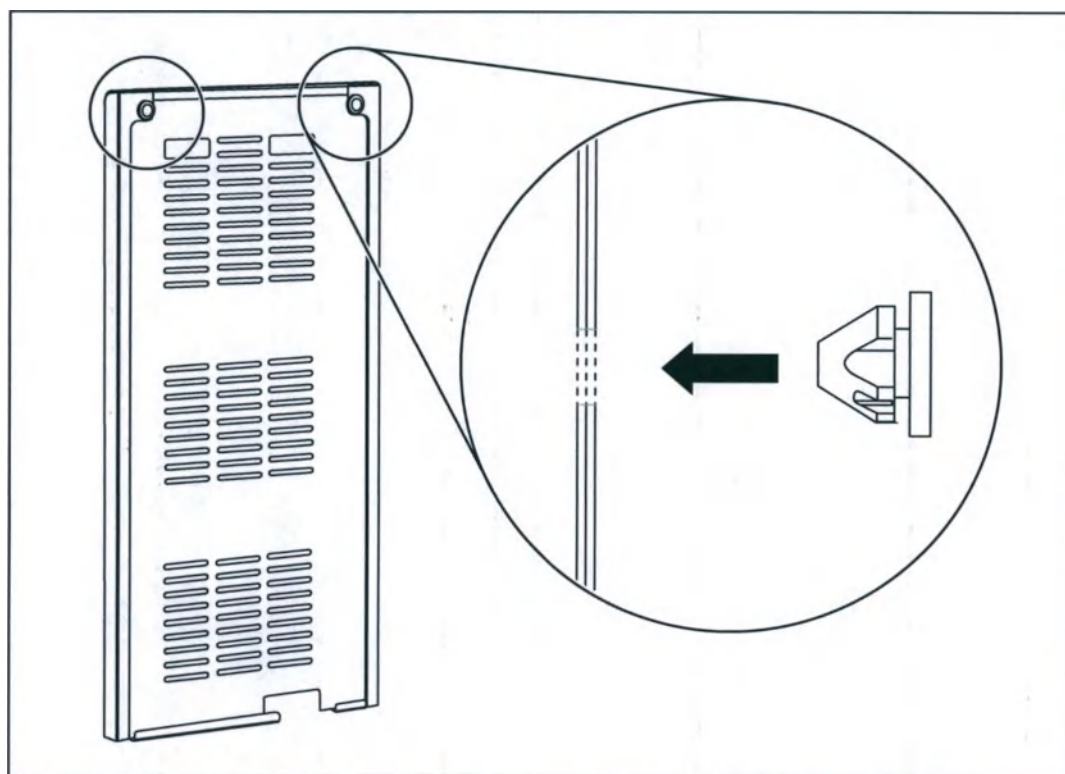2. Insert the new snaploc coupling into the hole and press firmly until retained.

12 Troubleshooting

CAUTION

Risk of injury

Should any irregularity be suspected, stop using the product immediately and refer to the table below. If the problem cannot be solved, repair is required, please contact the Endoscope Service Center, a service center designated by Olympus, an Olympus branch or a sales office.

Problem	Possible cause	Remedy
No power to separation transformer	Separation transformer is not connected to the mains electricity	Connect separation transformer to the mains electricity.
	No power at mains electricity	Check power supply to mains electricity.
Circuit breaker(s) 'trip'	Combined rating of medical electrical equipment connected to the outlet sockets exceeds the maximum VA rating of the separation transformer.	Check combined rating of medical electrical equipment is within the maximum limit, disconnect medical electrical equipment if necessary. Reset circuit breaker(s) by pressing in power switch / circuit breaker button.
	Operating the power button several times in rapid succession.	Reset circuit breaker(s) tip by pressing in power switch / circuit breaker button.
Electromagnetic interference.	Electromagnetic disturbance caused by the proximity of sensitive medical electrical equipment to an electrical field.	Relocate the sensitive medical electrical equipment to a higher shelf.
	Electromagnetic disturbance caused by damaged internals of the WM-NP3 Mobile Workstation.	Check the system for damaged cables or connections.
	Electromagnetic disturbance caused by unwanted signal transmission.	Investigate the source of the unwanted signal. Relocate the sensitive medical electrical equipment to a higher shelf.
No power to mounted equipment	Power cable connector is loose or not connected to the separation transformer	Check power cable connector.
	Power cable of the medical electrical equipment is not connected to mains electricity	Connect medical electrical equipment cables to the mains electricity.
Castor noise	Objects in braked castor	Check braked castors for damage or foreign objects.

13 End of life

Equipment which has reached end of its useful life must be decontaminated to local regulations before disposal and disposed of according to local regulations. Refer to "Disposal" on page 65 for additional information.

Refer to "Separation transformer" on page 78 for the disposal of the separation transformer.

14 Glossary

AP

Flammable anesthetic proof mixtures with air

APG

Flammable anesthetic proof mixtures with oxygen or nitrous oxide

EMC

Electromagnetic compatibility (EMC) is the ability of electrical equipment and systems to function satisfactorily in their electromagnetic environment without introducing intolerable electromagnetic disturbances to anything in that environment, such as electromagnetic interference (EMI) or even physical damage in operational equipment.

GI

Gastrointestinal

HF

High Frequency. Frequencies of alternating current above 100 kHz. High frequencies are used in electrosurgery, as lower frequencies provoke excessive neuromuscular stimulation and perhaps electrocution. High frequencies are also called radio frequencies (RF), as frequencies above 100 kHz can also transmit radio signals.

Leakage current

Leakage currents are current flows that pass another pathway than the initially intended electric circuit.

LCD

Liquid Crystal Display

Medical electrical equipment

Abbreviation for Medical electrical equipment. Electrical equipment having an applied part or transferring energy to or from the patient or detecting such energy transfer to or from the patient.

Reprocessing

Reprocessing is the process that is carried out on a used medical device in order to allow its safe reuse including cleaning, disinfection and, if applicable, sterilization and related procedures. Reprocessing also includes the testing and restoring of the technical and functional safety of the used device.

RF

Radiofrequency. Refer to the glossary entry "High frequency".

Separation transformer

Device to provide a separated power supply to electrical equipment to ensure system leakage currents are maintained within specified limits.

UDI

Unique Device Identifier. Abbreviation for Unique Device Identifier. The Unique Device Identifier is a label on the medical device with human- and machine-readable information, which allows to uniquely identify the device through distribution and use.

URO

Urology

WM-NP3

Workstation Mobile Narrow Platform 3

15 Index

B

Braked castors21

C

Cable management.....47

Castors.....21

Contact.....8

D

Disposal.....64

E

Electrical equipment.....43

Electromagnetic compatibility.....72

Equipment.....43

F

Fitting kit.....15

H

Height.....76

I

Irrigation Tube Holder.....27

K

Keyboard.....42

L

Load capacity.....76

M

Maneuvering.....51

Medical electrical equipment.....43

P

Potential equalization terminals.....	68
Power cable.....	20, 65
Protective earth terminal	66

R

Rear panel.....	29
Replacement	79
Reprocessing agent	57

S

Scope of delivery.....	14
Scope Pole	31
Separation transformer.....	44, 77
Shelves.....	23
Side handles	22
Spare parts.....	79
Specifications	75
Standard sets	11

T

Troubleshooting.....	98
Tubings.....	27

W

Warranty.....	20
Water bottle	29
Weight	76
Winding	67

OLYMPUS

OLYMPUS MEDICAL SYSTEMS CORP.

2851 Ishikawa-cho, Hachioji-shi, Tokyo 182-8507, Japan
Fax: (042) 646-2429 Telephone: (042) 642-2111

OLYMPUS EUROPA SE & CO. KG

Wendenstrasse 20, 20097 Hamburg, Germany
(Letters) Postfach 10 49 08, 20034 Hamburg, Germany
Fax: (040) 23773-4666 Telephone: (040) 23773-0

OLYMPUS AMERICA INC.

3500 Corporate Parkway, P.O. Box 610, Center Valley, PA
18034-0610, U.S.A.
Fax: (484) 896-7128 Telephone: (484) 896-5000

KEYMED (MEDICAL AND INDUSTRIAL EQUIPMENT) LTD.

KeyMed House, Stock Road, Southend-on-Sea, Essex SS2 5QH, United Kingdom
Fax: (01702) 465677 Telephone: (01702) 616333

OLYMPUS SINGAPORE PTE LTD.

438B, Alexandra Road, Alexandra Technopark Bldg B,
Fax: 6834-2436 Telephone: 6834-0010 #03-07712, Singapore 119968

OLYMPUS (BEIJING) SALES & SERVICE CO., LTD.

A8F, Ping An International Financial Center, No. 1-3, Xinyuan South Road,
Chaoyang District, Beijing, 100027 P.R.C.
Fax: (86)10-5876-1299 Telephone: (86)10-5819-9000

OLYMPUS MOSCOW LIMITED LIABILITY COMPANY

Elektronovodskaya str. 27 bld.8, 107023 Moscow, Russia
Fax: (7)495-826-7072 Telephone: (7)495-826-7077

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Mass, g	248 each
<i>Horizontal endoscope pole</i>	
Dimensions (HxDxW), mm	192.7 x 25.4 x 562.9
Mass, kg	1.1
<i>Vertical endoscope pole</i>	
Dimensions (HxDxW), mm	1111 x 25.4 x 214
Mass, kg	1.2
<i>Endoscope connector holder</i>	
Dimensions (HxDxW), mm	46 x 102 x 284
Mass, g	231
<i>Endoscope support</i>	
Dimensions (HxDxW), mm	125.8 x 68 x 68 each
Mass, g	35 each
<i>Scope pole mount</i>	
Dimensions (HxDxW), mm	197 x 50 x 236
Mass, g	785
<i>M6 T-nuts</i>	
Dimensions (LxWxH), mm	21.9 x 13.25 x 7.6
Mass, g	10g
<i>M6 button-head screws</i>	
Dimensions (LxWxH), mm	10.35 Dia x 23.4
Mass, g	5
<i>M6 countersunk screws</i>	
Dimensions (LxWxH), mm	11.5 Dia x 16.2
Mass, g	4
<i>M6 crinkle washer</i>	
Dimensions (LxWxH), mm	11.3 Dia x 0.7
Mass, g	<1
<i>M6 plain washer</i>	
Dimensions (LxWxH), mm	11.8 Dia x 1.6
Mass, g	1g

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Tolerance is $\pm 5\%$ unless otherwise specified.

Overall dimensions of Accessories

Characteristics	Value
EUS connector holder (MAJ-1635)	
Dimensions (LxWxH), mm	137.5 x 361 x 155
Mass, kg	1.1
Single CO₂ holder (up to Ø140 mm) (MAJ-1639)	

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<i>Fitting kit</i>	
<i>M6 crinkle washer</i>	
Dimensions (HxDxW), mm	11.3 Dia x 0.7
Mass, g	<1
<i>M6 plain washer</i>	
Dimensions (HxDxW), mm	11.8 Dia x 1.6
Mass, g	1
<i>M6 button-head screw</i>	
Dimensions (HxDxW), mm	10.35 Dia x 23.4
Mass, g	5
<i>4 mm hex key</i>	
Dimensions (HxDxW), mm	141 x 28 x 4.5
Mass, g	17
<i>Sterile Water Holder (MAJ-2166)</i>	
Dimensions (LxWxH), mm	99 x 110 x 101.5
Mass, g	310
<i>Nurses Control Arm (MAJ-2167)</i>	
Dimensions (HxDxW), mm	375 x 110 x 160
Mass, kg	1,7
Maximum monitor weight	9 kg
Maximum monitor size (HxDxW), mm	420 x 131 x 586
<i>Fitting kit</i>	
<i>M6 T-nuts</i>	
Dimensions (HxDxW), mm	21.9 x 13.25 x 7.6
Mass, g	10
<i>M6 flat-head screws</i>	
Dimensions (HxDxW), mm	11.8 Dia x 25
Mass, g	5
<i>4 mm hex key</i>	
Dimensions (HxDxW), mm	141 x 28 x 4.5

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Mass, g	17
<i>Monitor fixing screws</i>	
Dimensions (HxDxW), mm	7.8 Dia x 13
Mass, g	2
Side Shelf (MAJ-2173)	
Dimensions (HxDxW), mm	508 x 216 x 365
Mass, g	860
Maximum weight capacity	5 kg
<i>Shelf Strap</i>	
Dimensions (HxDxW), mm	560 x 20 x 32
Mass, g	29
Fitting kit	
<i>M6 button-head screws</i>	
Dimensions (HxDxW), mm	10.35 Dia x 23.4
Mass, g	5
<i>M6 crinkle</i>	
Dimensions (HxDxW), mm	11.3 Dia x 0.7
Mass, g	<1
<i>M6 plain washers</i>	
Dimensions (HxDxW), mm	11.8 Dia x 1.6
Mass, g	1
<i>4 mm hex key</i>	
Dimensions (HxDxW), mm	141 x 28 x 4.5
Mass, g	17
LCD Monitor Arm 12-14 kg (MAJ-2217)	
Dimensions (HxDxW), mm	380 x 200 x 482
Mass, kg	7.1
Maximum monitor weight	12 kg - 14 kg
Fitting kit	
<i>M8 T-nuts</i>	
Dimensions (HxDxW), mm	21.9 x 13.25 x 7.6

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Mass, g	10
<i>M8 screws</i>	
Dimensions (HxDxW), mm	13.1 Dia x 28
Mass, g	15
<i>M8 crinkle washers</i>	
Dimensions (HxDxW), mm	14.7 Dia x 1.2
Mass, g	<1
<i>6 mm hex key</i>	
Dimensions (HxDxW), mm	140 x 38 x 6.76
Mass, g	39
<i>Monitor fixing screws</i>	
Dimensions (HxDxW), mm	7.8 Dia x 13
Mass, g	2

Tolerance is $\pm 5\%$ unless otherwise specified.

Key accessories, medical devices or non-medical devices, to be used in combination with medical device submitted for registration

The WM-NP3 will be compatible with the following capital equipment configurations:

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Other Equipment:

Name	Registration certificate number
IMH-10 Image Management Hub	FSZ 2011/10652
IMH-20 Image Management Hub	
nCare	-
ESG-100 HF Electrosurgical Generator	FSZ 2009/03927
ESG-150 Electrosurgical Generator	-
ESG-200 Electrosurgical Generator	RZN 2015/3465
Sony UP25-MD Video Colour Printer	-

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- 1.8.2 Side Handle (right) – 1 pc
- 1.8.3 Fitting kit – 1 pc, containing:
 - 1.8.3.1 M6 T-nuts – 4 pc
 - 1.8.3.2 M6 countersunk screws – 4 pc
 - 1.8.3.3 4 mm hex key – 1 pc
- 1.9 IEC lead set (optional):
 - 1.10.1 Cable universal IEC 320 0.8 m – 3 pc
 - 1.10.2 Cable universal IEC 320 1.5 m – 5 pc
 - 1.10.3 Cable universal IEC 320 M-FEM 2.5 m – 1 pc
- 1.10 Water bottle bracket
- 1.11 Instructions for use
- 2. WM-NP3 Mobile Workstation Basic:
 - 2.1 WM-NP3 Mobile Workstation
 - 2.2 Workstation fitting kit:
 - 2.2.1 4 mm hex key – 1 pc
 - 2.2.2 6 mm hex key – 1 pc
 - 2.2.3 Cable ties – 12 pc
 - 2.3 Separation transformer (MAJ-2157)
 - 2.4 Irrigation Tube Holder (MAJ-2160) (optional)
 - 2.5 IEC lead set:
 - 2.5.1 Cable universal IEC 320 0.8 m – 3 pc
 - 2.5.2 Cable universal IEC 320 1.5 m – 5 pc
 - 2.5.3 Cable universal IEC 320 M-FEM 2.5 m – 1 pc
 - 2.6 Water bottle bracket
 - 2.7 Instructions for use
- 3. WM-NP3 Mobile Workstation URO:
 - 3.1 WM-NP3 Mobile Workstation
 - 3.2 Workstation fitting kit:
 - 3.2.1 4 mm hex key – 1 pc
 - 3.2.2 6 mm hex key – 1 pc
 - 3.2.3 Cable ties – 12 pc
 - 3.3 Separation transformer (MAJ-2157)

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

	Mains power ON / OFF switch
	MDR EU regulatory requirement (Class I)
	Do not push the workstation when the brakes are applied
	Max laden weight
	Unit complies with China RoHS 2 requirements
	Stability warning
	Unit complies with METI-PSE regulations – only applicable for Japan
	KC marking indicating that the device complies with the requirements of the Law Radio Frequencies of Korea (Paragraph 3, Article 58-2)
	CE marking indicating that the device is in conformity with the applicable requirements set out in applicable European Union harmonization legislation.
	Consult instructions for use
	eIfU - Instructions for use available online
	Catalogue number

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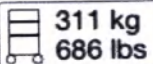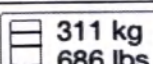	Input
	Output
	Alternating current
	Weight
	Remote switch
	Non-inherently short circuit proof separation transformer
	According to European Directive 2012/19/EU on Waste Electrical and Electronic Equipment, this symbol indicates that the product must not be disposed off as unsorted municipal waste, but should be collected separately.
	Indicates the unit has been independently assessed for compliance with UL and CSA standards as a Recognized Component of The WM-NP3 Mobile Workstation (MAJ-2156 only).
	European Authorized Representative
	Translation
	Importer
	UKCA marking indicating that the device is in conformity with the applicable requirements set out in

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© KeyMed 2021
(Keymed (Medical & Industrial Equipment) Ltd)