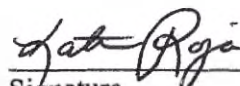
Medtronic

18000 Devonshire Street
Northridge, CA 91325 USA
www.medtronic.com
Tel 818-576-5211

ATTESTATION STATEMENT

As a responsible member of Medtronic MiniMed, I hereby attest that I believe the following Appendix to IFU attached herein, to be a true and accurate copy.

Katherine Rojas
Regulatory Affairs Specialist

Signature

A notary public or other officer completing this certificate verifies only the identity of the individual who signed the document to which this certificate is attached, and not the truthfulness, accuracy, or validity of that document

State of California

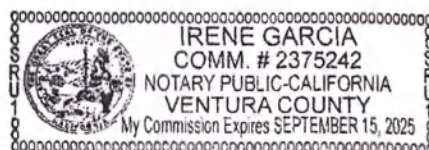
County of VENTURA

On JUNE 21, 2022 before me, IRENE GARCIA, NOTARY PUBLIC,
(date) (insert name and title of the officer)

personally appeared KATHERINE VIRGINIA ROJAS, who proved to me on the basis of satisfactory evidence to be the person(s) whose name(s) is/are subscribed to the within instrument and acknowledged to me that he/she/they executed the same in his/her/their authorized capacity(ies), and that by his/her/their signature(s) on the instrument the person(s), or the entity upon behalf of which the person(s) acted, executed the instrument. I certify under PENALTY OF PERJURY under the laws of the State of California that the foregoing paragraph is true and correct.

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Transmitter Guardian Connect (MMT-7821, MMT-7821L) is part of the Guardian Connect Kit (MMT-7820, MMT-7820L). The technical characteristics of the MMT-7821 model are identical in all respects to the model of the transmitter included in the MMT-7820 kit.

Transmitter (MMT-7821) specifications	
Protection class	IP48
Type of applied part	BF
Software Version	V2.1B
Software safety class	B
Electrical shock protection class	Class II
Software release date	14 Oct 2016
Dimensions	3.5 cm (w) x 4.1 (l) x 1.3 cm (h)
Weight	4.75 g
Sterilization	Non-sterile device
Shelf Life	Six (6) months
Operational Life	One (1) year
Battery Life	Seven (7) days following a full charge
Battery	Lithium Polymer
Radiofrequency/Rf Telemetry	BLE v4.0 wireless (2.4 GHz)
Compatible devices	
CSS7200, CSS7201	Guardian Connect Application
MMT-7008	Glucose Sensor (Enlite) (license no. RZN 2013/785 от 04.07.2013)
Transmitter (MMT-7821L) Specifications	
Protection class	IP48
Type of applied part	BF

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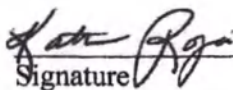
ATTESTATION STATEMENT

As a responsible member of Medtronic MiniMed, I hereby attest that I believe the following Instructions for Use attached herein, to be a true and accurate copies.

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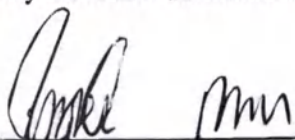
State of California

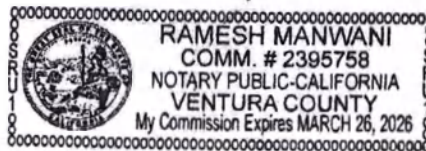
County of Ventura

On 5 July, 2022 before me, Ramesh Manwani (Notary Public)
(date) (insert name and title of the officer)

personally appeared Katherine Rojas, who proved to me on the basis of satisfactory evidence to be the person(s) whose name(s) is/~~are~~ subscribed to the within instrument and acknowledged to me that ~~he~~/she/~~they~~ executed the same in ~~his~~/her/~~their~~ authorized capacity(ies), and that by ~~his~~/her/~~their~~ signature(s) on the instrument the person(s), or the entity upon behalf of which the person(s) acted, executed the instrument. I certify under PENALTY OF PERJURY under the laws of the State of California that the foregoing paragraph is true and correct.

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Medtronic

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Africa: Medtronic South Africa and Southern Africa

Office Reception Tel: +27(0) 11 260 9300
Diabetes: 24/7 Helpline: 0800 633 7867
Sub-Sahara 24/7 Helpline: +27(0) 11 260 9490

Albania: Net Electronics Albania
Tel: +355 697070121

Argentina: Corpomedica S.A.
Tel: +(11) 4 814 1333
Medtronic Directo 24/7:
+0800 333 0752

Armenia: Exiol LLC
Tel: +374 98 92 00 11
or +374 94 38 38 52

Australia: Medtronic Australasia Pty. Ltd.
Tel: 1800 668 670

Bangladesh: Sonargaon Healthcare Pvt Ltd.
Mobile: (+91)-9903995417
or (+880)-1714217131

Belarus: Zarga Medica
Tel: +375 29 625 07 77
or +375 44 733 30 99
Helpline: +74995830400

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Malaysia: Medtronic International Ltd.
Tel: +603 7946 9000

Middle East and North Africa: Regional Office
Tel: +961-1-370 670

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Medtronic Australasia Pty. Ltd.
Tel.: 1800 668 670

Bangladesh:

Sonargaon Healthcare Pvt Ltd.
Mobile: (+91)-9903995417
or (+880)-1714217131

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USA:

Medtronic Diabetes Global
Headquarters
24-Hour Technical Support:
+1-800-646-4633
To order supplies: +1-800-843-6687

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Introduction

The tester (MMT-7736) is an accessory item included with select Medtronic Diabetes continuous glucose transmitters.

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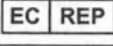
Disposal

Dispose of the tester according to local regulations for electronic waste.

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Icon table

	Consult instructions for use
	Caution
	Manufacturer
	Date of manufacture
(2x)	Two per container/package
	Catalogue number
	Batch code
	Storage temperature limit
	Storage humidity limit
	Importer
	Do not dispose of this product in unsorted municipal waste stream
	Medical device
	Authorized representative in the European Community
	Signifies European conformity.

© 2019 Medtronic MiniMed, Inc. All rights reserved. MiniMed™ is a trademark of Medtronic MiniMed, Inc.

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Introduction

The tester (MMT-7736L) is an accessory item included with select Medtronic Diabetes continuous glucose transmitters.

Indications for use

The tester is used to test the transmitter. The tester is also used to form a waterproof seal that protects the transmitter connector while the transmitter is being cleaned.

Contraindications

None known.

Clinical benefits

Not applicable. The tester is an accessory that tests the transmitter and protects it during cleaning.

User safety

Warnings

Keep the tester out of reach of children. The tester may pose a choking hazard that can result in serious injury or death.

Do not use the tester if it comes into contact with blood. Contact with blood may cause infection. If the tester comes into contact with blood, dispose directly into a sharps container.

Precautions

Read the entire transmitter user guide to prevent misuse of the tester. Using the tester before reading the transmitter user guide can lead to misuse of the transmitter.

Risks and side effects

The tester may pose a choking hazard to children, which can result in serious injury or death.

Hazardous substances

None.

Allergens

Internal components contain nickel.

Using the tester

Refer to the transmitter user guide for instructions.

Maintenance

Cleaning

Not applicable.

Storage

Store the tester at room temperature.

Disposal

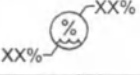
Dispose of the tester according to local regulations for electronic waste disposal.

Technical specifications

The tester has an expected lifetime of one year when connected to the transmitter once per sensor replacement period. If physical damage or discoloration is observed, the tester should be replaced.

Approximate dimensions	Approximate weight
3.8 x 2.0 x 1.0 centimeters (1.5 x 0.8 x 0.3 inches)	3.0 grams (0.11 ounces)

Icon table

Symbol	Meaning
	Manufacturer
	Authorized representative in the European Community.
	Catalogue number
	Storage temperature limits
	Consult instructions for use
	Medical device
	Importer
	Date of manufacture
(2x)	Two per container/package
	Batch code
	Storage humidity limits
	Do not dispose of this product in unsorted municipal waste stream
	Caution
	Conformité Européenne (European Conformity). This symbol means that the device fully complies with applicable European Union Acts.

© 2019 Medtronic. All rights reserved. Medtronic, Medtronic logo, and Further, Together are trademarks of Medtronic. All other brands are trademarks of a Medtronic company.

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USA:

Medtronic Diabetes Global
Headquarters
24-Hour Technical Support:
+1-800-646-4633
To order supplies:
+1-800-843-6687

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Introduction

The charger (MMT-7715) is an accessory intended for use with Medtronic continuous glucose monitoring (CGM) transmitter devices. Transmitters contain a non-replaceable, rechargeable battery. The charger is used to recharge a transmitter's internal battery.

Indications for use

The charger is indicated for use with Medtronic CGM transmitters. Refer to the system user guide for additional information, such as the intended age group.

Contraindications

No contraindications are associated with charger use. For contraindications related to continuous glucose monitoring, see the system user guide.

Clinical benefits

Not applicable. The charger is an accessory that is used to recharge the internal battery of Medtronic transmitters.

User safety

Warnings

Keep the charger out of reach of children. The charger contains small parts and may pose a choking hazard that can result in serious injury or death.

For questions or concerns related to charger use, contact a local Medtronic support representative for assistance.

For medical questions or concerns, or to report a serious incident, contact a healthcare professional.

Precautions

Consult the system user guide for precautions related to continuous glucose monitoring.

Risks and side effects

The charger contains small parts and may pose a choking hazard that can result in serious injury or death.

Hazardous substances

None.

Allergens

None known.

Using the charger

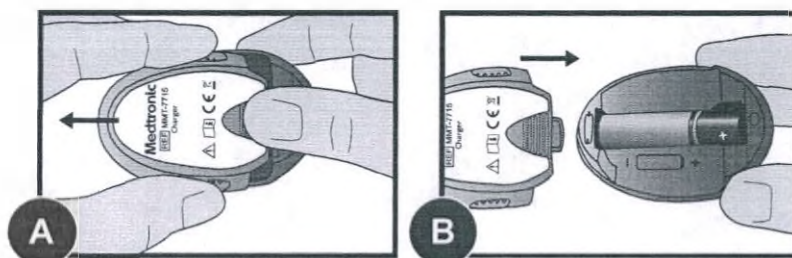
The charger has a green light that shows the charging status and a red light that communicates any problems during charging. If you see a red light, see the **Troubleshooting** section. The charger requires one AAA alkaline battery.

Note: If the battery is installed incorrectly, or is low, the charger will not work. Repeat the battery installation steps using a new battery.

Installing a battery in the charger

To install a battery in the charger:

1. Push the battery cover in and slide it off (as shown in image A).
2. Insert a new AAA alkaline battery. Make sure the + and - symbols on the battery align with these same symbols shown on the charger.
3. Slide the cover back on the charger until it clicks into place (as shown in image B).

Troubleshooting

Problem	Likely Cause(s)	Resolution
A device is connected to the charger and no lights come on.	The charger battery has no power or no battery is inserted.	If there is no damage to the connector pins, replace the battery in the charger.
During charging, the flashing green light on the charger turns off and a longer flashing red light appears on the charger.	The charger battery is low on power.	Replace the battery in the charger.
During charging, a mix of quick and long flashing red lights appear on the charger.	The charger <i>and</i> the device are low on power.	Replace the battery in the charger and fully charge the device.

Maintenance

Cleaning

CAUTION: Do not immerse the charger in water or any other cleaning agent. The charger is not waterproof. Water can damage the charger and cause the device to malfunction.

To clean the charger:

1. Wash your hands thoroughly.
2. Use a damp cloth with mild cleaning solution, such as a dish washing detergent, to clean any dirt or foreign material from the outside of the charger. Never use organic solvents such as paint thinner or acetone to clean the charger.
3. Place the charger on a clean, dry cloth and air dry for two to three minutes.

Storage

Store the charger in a clean, dry location at room temperature.

Disposal

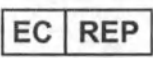
Dispose of the charger according to local regulations for electronic waste disposal.

Technical specifications

Operating Conditions	Temperature: 10° C to 40° C (50° F to 104° F) Relative Humidity: 30 % to 75 %, non-condensing
Storage Conditions	Temperature: -10° C to 50° C (14° F to 122° F) Relative Humidity: 10 % to 95 %, non-condensing

Icon Table

	Consult instructions for use
	Caution: consult instructions for use for important warnings or precautions not found on the label
	Batch code
	Catalogue number
	One per container/package

	Storage humidity limits
	Storage temperature limits
	Manufacturer
	Date of manufacture
	Medical device
	Importer
	Do not dispose of this product in unsorted municipal waste stream
	Authorized representative in the European Community
	Conformité Européenne (European Conformity). This symbol means that the device fully complies with applicable European Union Acts

© 2019 Medtronic. Medtronic, Medtronic with rising man logo, and Medtronic logo are trademarks of Medtronic.

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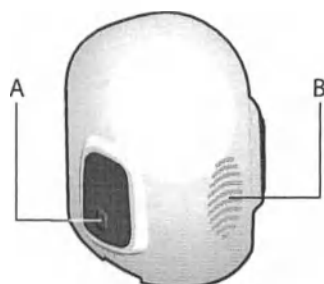
Introduction

The one-pressserter (MMT-7512) is used as an aid to insert Medtronic MiniMed™ glucose sensors.

- A. Bump (on both buttons)
- B. Thumbprint marking

Indications for use

The sserter is indicated for single-patient use. It is not intended for use by multiple patients. Refer to the system user guide for additional information, such as the intended age group.

Contraindications

No contraindications are associated with sserter use. For contraindications related to continuous glucose monitoring, see the system user guide.

Clinical benefits

Not applicable. The sserter is an accessory that helps to insert compatible Medtronic sensors.

User Safety

Warnings

Read this entire user guide prior to sensor insertion. The sserter does not work the same way as other Medtronic insertion devices. Failure to follow directions may result in improper insertion, pain, or injury.

For additional warnings related to sensor use, see the sensor user guide.

Keep the sserter out of reach of children. The sserter contains small parts and may pose a choking hazard that can result in serious injury or death.

Do not make therapy decisions based on sensor glucose values. Sensor glucose and blood glucose values may differ. If a sensor glucose reading is low or high, or symptoms of low or high blood glucose appear, use a BG meter to confirm blood glucose prior to making therapy decisions. Failure to check blood glucose levels with a BG meter prior to making therapy decisions can result in the infusion of too much or too little insulin, which can cause hypoglycemia or hyperglycemia.

Healthcare professionals and caregivers:

- A retractable needle is attached to the sensor. Minimal bleeding may occur. Always wear gloves to insert the sensor
- Cover the sensor with sterile gauze to remove the needle housing from the sensor

Place the needle housing directly into a sharps container after sensor insertion to prevent accidental needle stick injury.

Always inspect the sensor packaging for damage before use. Sensors are sterile and non-pyrogenic unless the package has been opened or damaged. If the sensor packaging is open or damaged, discard the sensor directly into a medical waste container. Use of a non-sterile sensor can result in infection at the insertion site.

For questions or concerns related to sensor use, contact a local Medtronic support representative for assistance.
For medical questions or concerns, contact a healthcare provider.

Precautions

Wash hands with soap and water prior to sensor insertion to reduce the risk of infection.

Inspect the insertion site for bleeding under, around, or on top of the sensor. If bleeding occurs, use a sterile gauze pad to apply steady pressure to the sensor for three minutes.

WARNING: Do not use non-sterile gauze to apply pressure to the sensor. Use of non-sterile gauze can cause infection at the insertion site.

If bleeding stops after three minutes of steady pressure, connect the transmitter to the sensor.

CAUTION: Do not connect the transmitter to the sensor if bleeding at the insertion site does not stop. Active bleeding can damage the transmitter and result in a loss of communication between the transmitter and the receiving device.

If bleeding does not stop after three minutes of steady pressure, remove the sensor, and continue to apply pressure until bleeding stops. Discard the used sensor directly in a sharps container after removal. Insert a new sensor at a different insertion site.

Note: Always use a new insertion site to insert a new sensor. Reuse of the same insertion site can cause redness, irritation, pain, bleeding, tenderness, and inflammation.

Do not clean, resterilize, or attempt to remove the needle from the needle housing. Needle stick injury may occur.

Do not reuse the sensor. Reuse of the sensor can cause damage to the sensor surface and lead to inaccurate sensor glucose values, irritation, and infection at the insertion site.

Rotate the sensor insertion site so that sites do not become overused.

Confirm that the sensor is firmly attached to the skin to prevent it from coming out of the insertion site. Factors that increase the likelihood of the sensor coming out include increased physical activity (particularly in younger patients) and improper taping technique.

Risks and side effects

The sensor contains small parts and may pose a choking hazard that can result in serious injury or death.

Side effects include discomfort and skin irritation at the insertion site.

Hazardous substances

None.

Allergens

None known.

Where to insert the sensor

Choose an insertion site with a sufficient layer of subcutaneous fat.

Note: Do not insert the sensor within 5.0 cm (2 inches) of the navel area to avoid discomfort at the insertion site and incomplete sensor adhesion.

For optimal sensor performance, and to prevent accidental sensor removal:

- Do not insert the sensor into muscle, tough skin, or scar tissue
- Avoid areas that are constrained by clothing or accessories
- Avoid areas that are located under a belt or waistband

- Avoid areas subjected to vigorous movement during exercise

For a diagram of recommended sensor insertion areas, see the sensor user guide.

Inserting the sensor

For instructions on how to insert the sensor, see the sensor user guide.

Maintenance

Cleaning

Clean the serter after every use or whenever debris, blood, or other contaminants are visible.

To clean the serter:

Note: Press and release the mechanism on the underside of the serter to fully rinse the device.

1. Rinse the serter under running tap water.
2. Submerge the serter in a solution of mild soap and warm water.
3. Use a soft-bristled brush to clean the serter until it is free of visible dirt.
4. Rinse the serter under running tap water until all visible soap is gone. Make sure to fully rinse all hard-to-reach areas.
5. Shake off any excess water. Place the serter upright on a clean cloth to air dry.

Storage

Store the serter in the released position at room temperature.

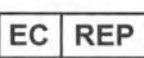
Disposal

Dispose of the serter according to local regulations.

Technical specifications

Approximate Dimensions	Approximate Weight
7.85 x 6.91 x 5.77 centimeters (3.09 x 2.72 x 2.27 inches)	112.5 grams (3.97 ounces)

Icon table

	Consult instructions for use
	Caution
	Batch code
	Catalogue number
	One per container/package
	Storage humidity limit
	Storage temperature limit
	Manufacturer
	Date of manufacture
	Importer
	Signifies European technical conformity
	Medical device
	Authorized representative in the European Community

©2019, Medtronic MiniMed, Inc. All rights reserved. MiniMed™ is a trademark of Medtronic MiniMed, Inc.

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United States

FCC ID: OH27821

This device complies with part 15 of the FCC Rules. Operation is subject to two conditions:

- (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Note: This device has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation.

This equipment generates, uses, and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this device does cause harmful interference to radio or television reception, which can be determined by turning the device off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the device and the receiver.
- Decrease the distance between the transmitter and the compatible mobile device to 6.1 meters (20 feet) or less.
- Increase the separation between the transmitter and the equipment that is receiving/interfering.

IMPORTANT: Do not change or modify the internal RF transmitter or antenna unless expressly approved by Medtronic Diabetes. Doing so could interfere with your ability to operate the device.

Canada

IC:3408B-7821

This device complies with Industry Canada license-exempt RSS standard(s). Operation is subject to the following two conditions:

- (1) this device may not cause interference, and (2) this device must accept any interference, including interference that may cause undesired operation of the device.

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Australia and New Zealand

This product conforms to Australia and New Zealand Radio Requirements.

Japan

This product conforms to Japan Radio Law Requirements.

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South Africa

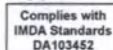
TA-2016/2275

This product conforms to South Africa Radio Requirements.

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Singapore

This product complies with IMDA Standards.

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(European Conformity)

EN Hereby, Medtronic declares that this Radio Equipment is in compliance with Directive 2014/53/EU.

The full text of the EU declaration of conformity is available at the following internet address:

www.MedtronicCompliance.eu

This product operates in the 2.4-2.4835 GHz band with a radiated output power of less than 1.57 mW e.i.r.p.

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© 2019 Medtronic MiniMed, Inc. All rights reserved. MiniMed™는 Medtronic MiniMed, Inc.의 상표입니다.
©2019 Medtronic MiniMed, Inc. All rights reserved. ミニメド™はMedtronic MiniMed, Inc.の商標です。
© 2019 Medtronic. All rights reserved. Medtronic, Medtronic 로고 및 Further, Together는 Medtronic의 상표입니다. 기타 모든 브랜드는 Medtronic 기업의 상표입니다.
©2019, Medtronic MiniMed, Inc. All rights reserved. MiniMed™는 Medtronic MiniMed, Inc.의 상표입니다.