

ПОДТВЕРЖДАЮ
Генеральный директор
АО «Трейдомед Инвест»



Соболева Э.Ж.
«15» ноября 2022 г.

**Фотографические изображения
на медицинское изделие**

«Система офтальмологическая электрофизиологическая, вариант исполнения Espion Profile»


Общий вид системы Espion Profile в конфигурации Gold



Тележка, номер модели D356, в составе:




- Столешница

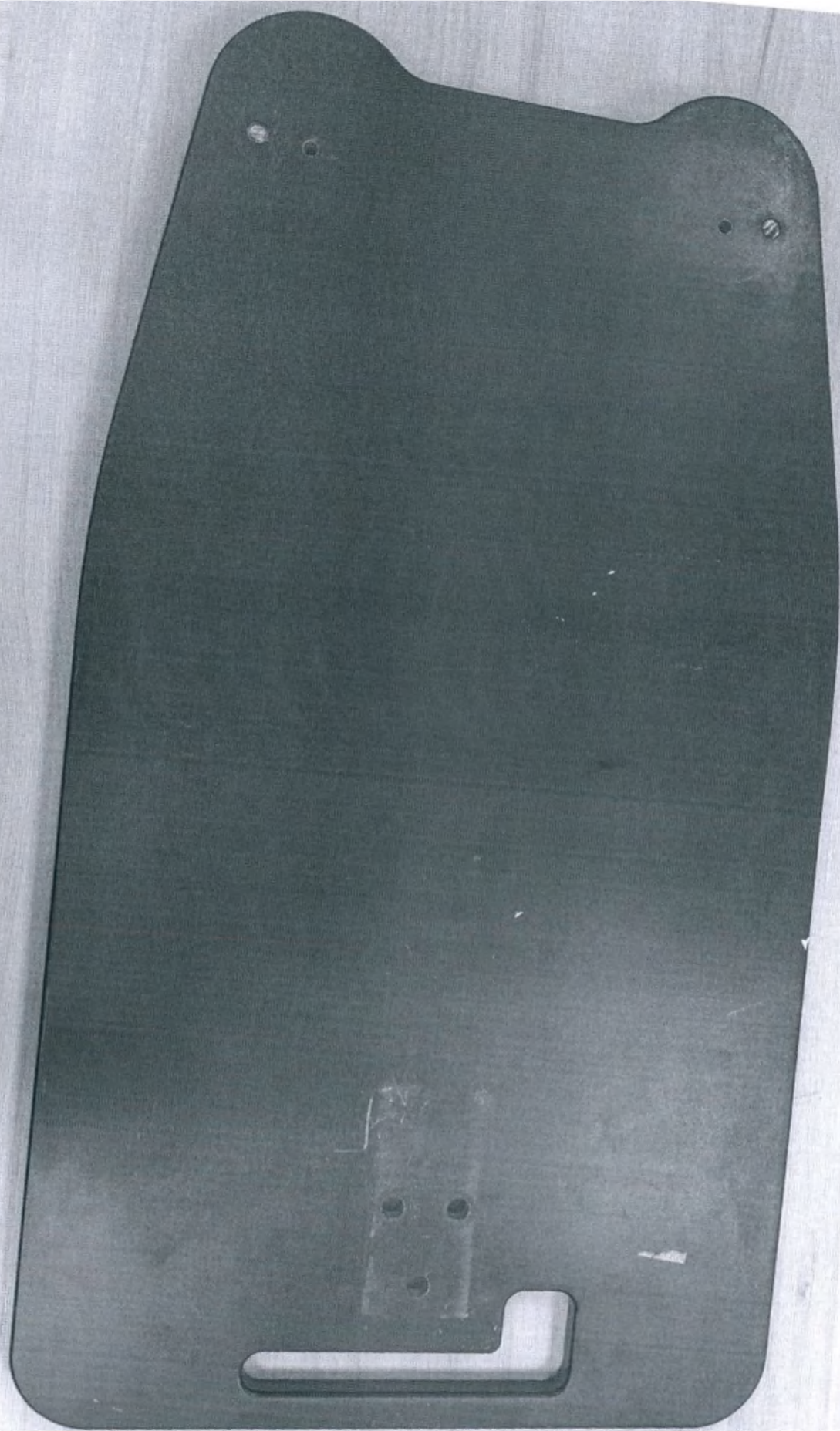





- Колесная рама – 2 шт.




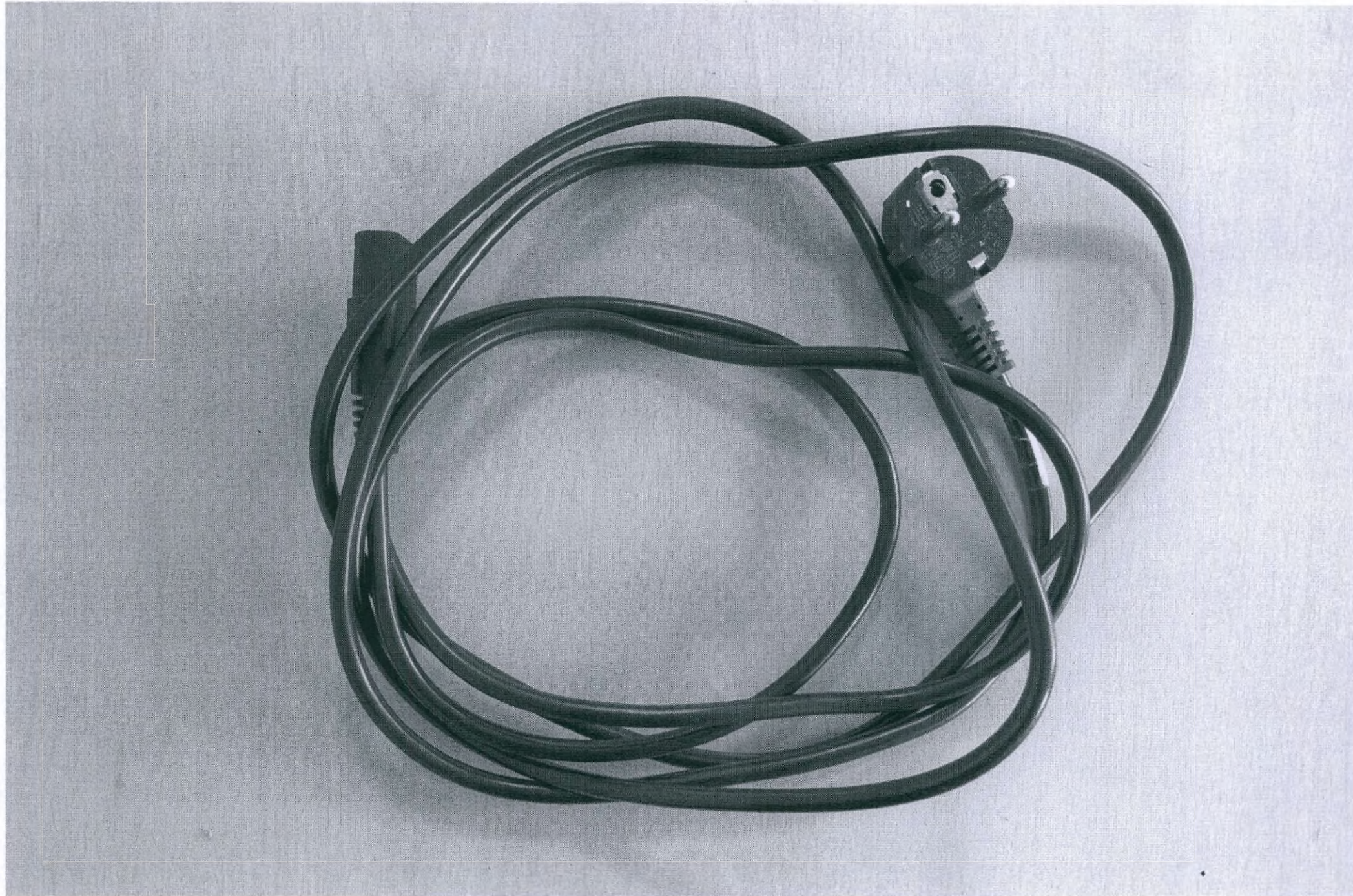

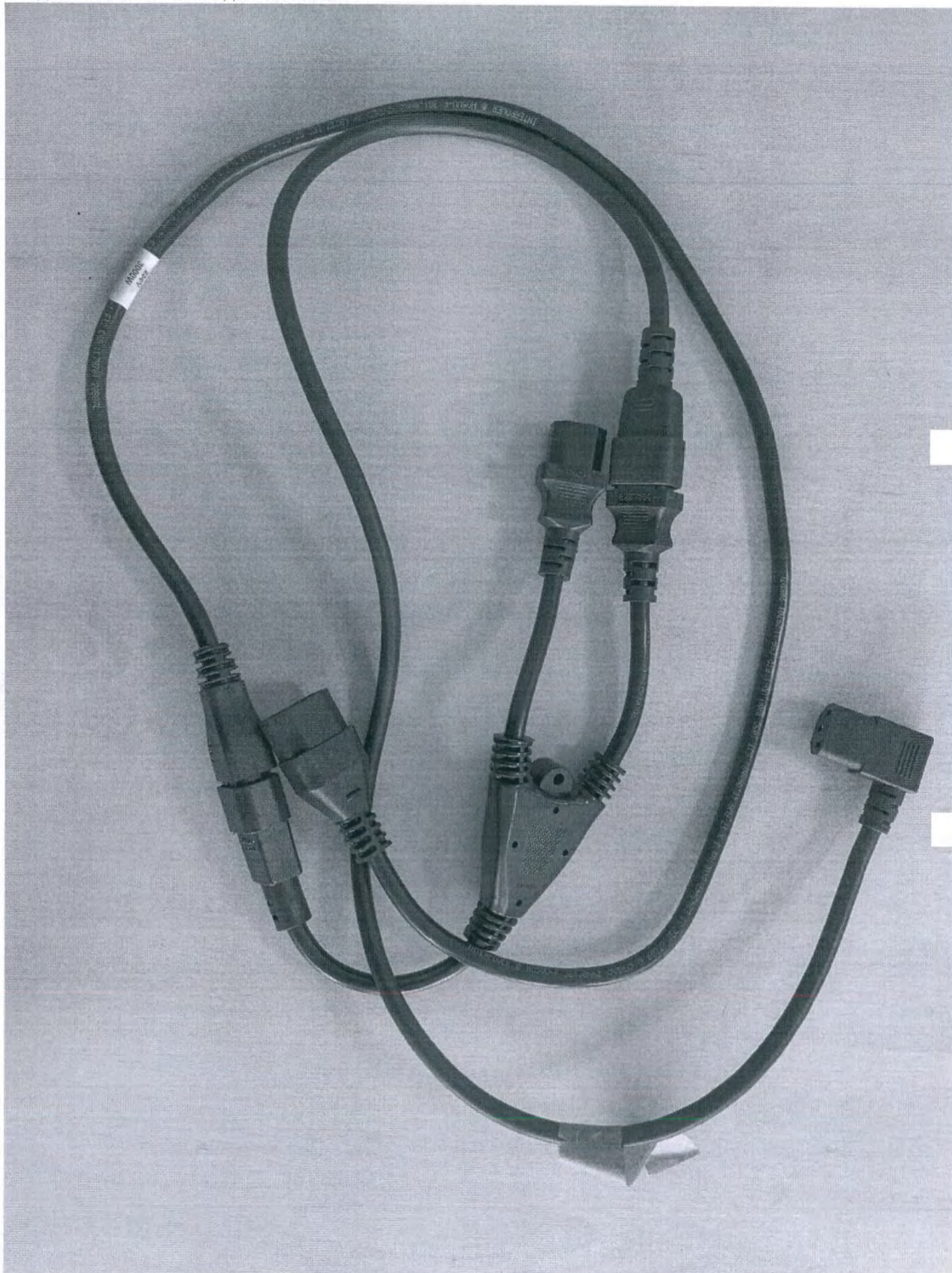



Кронштейн-столик для клавиатуры и мыши

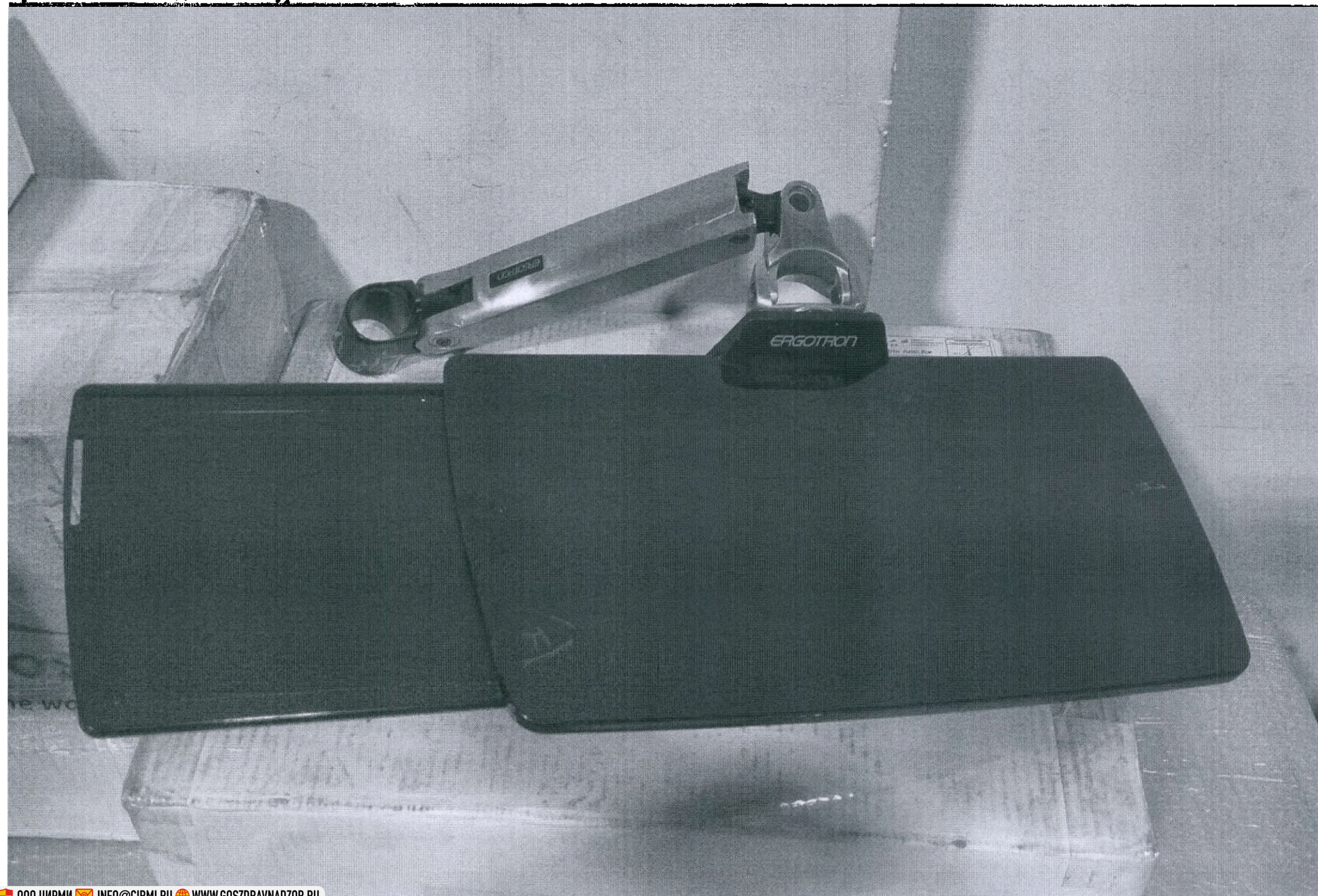


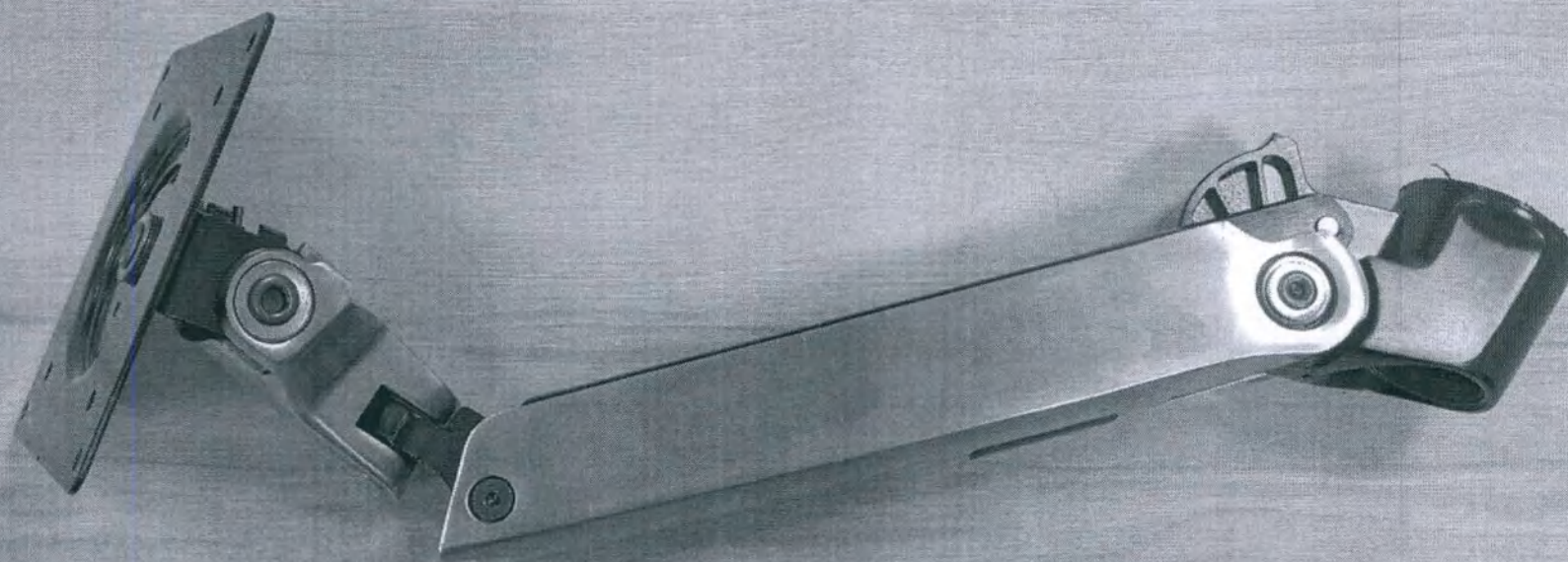



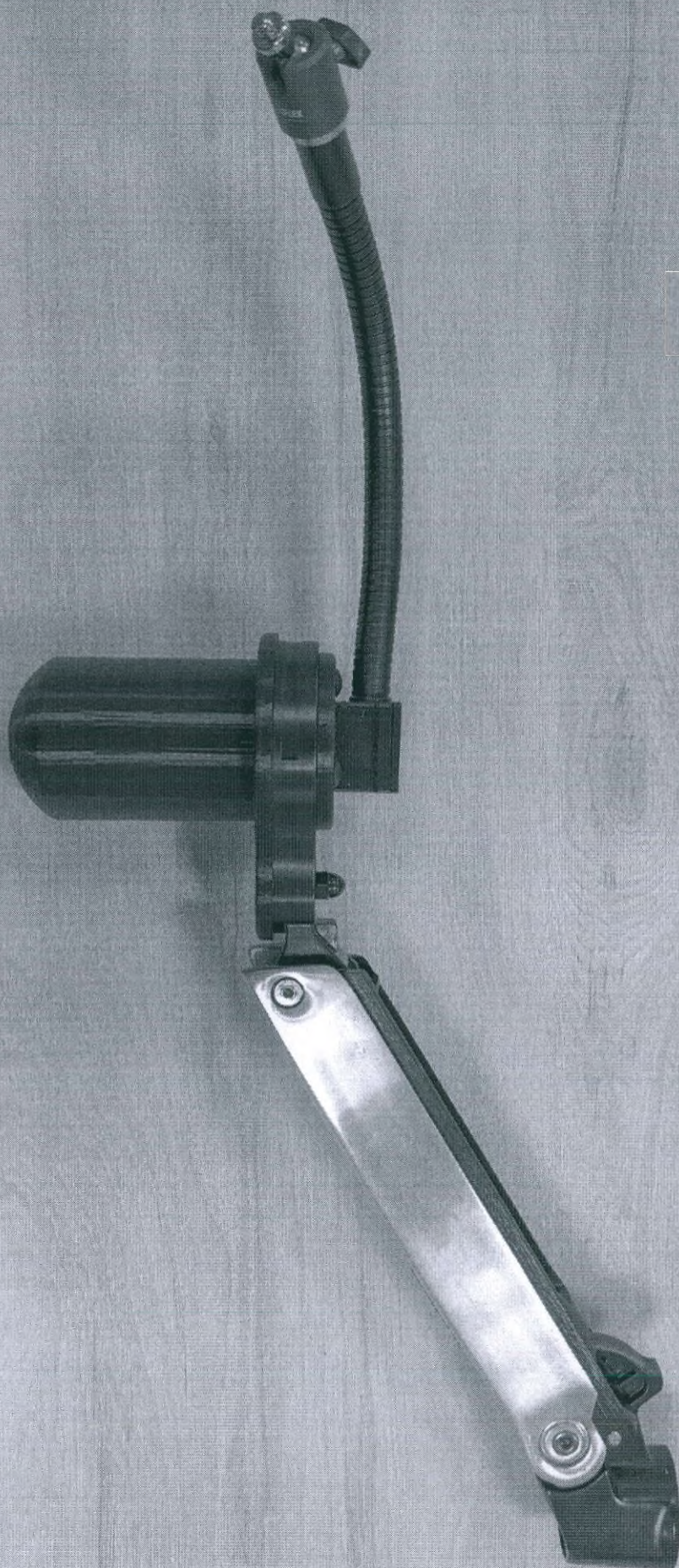


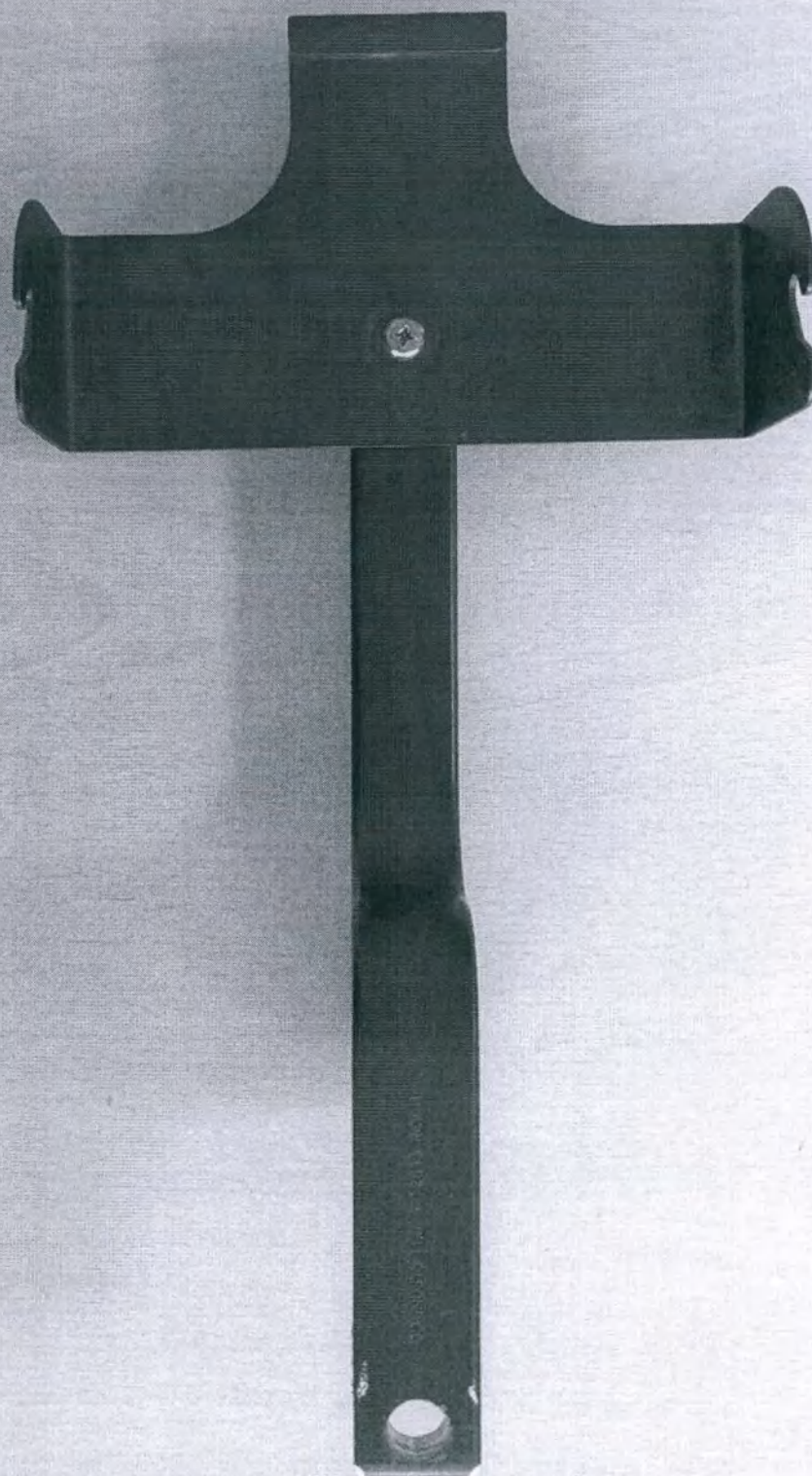


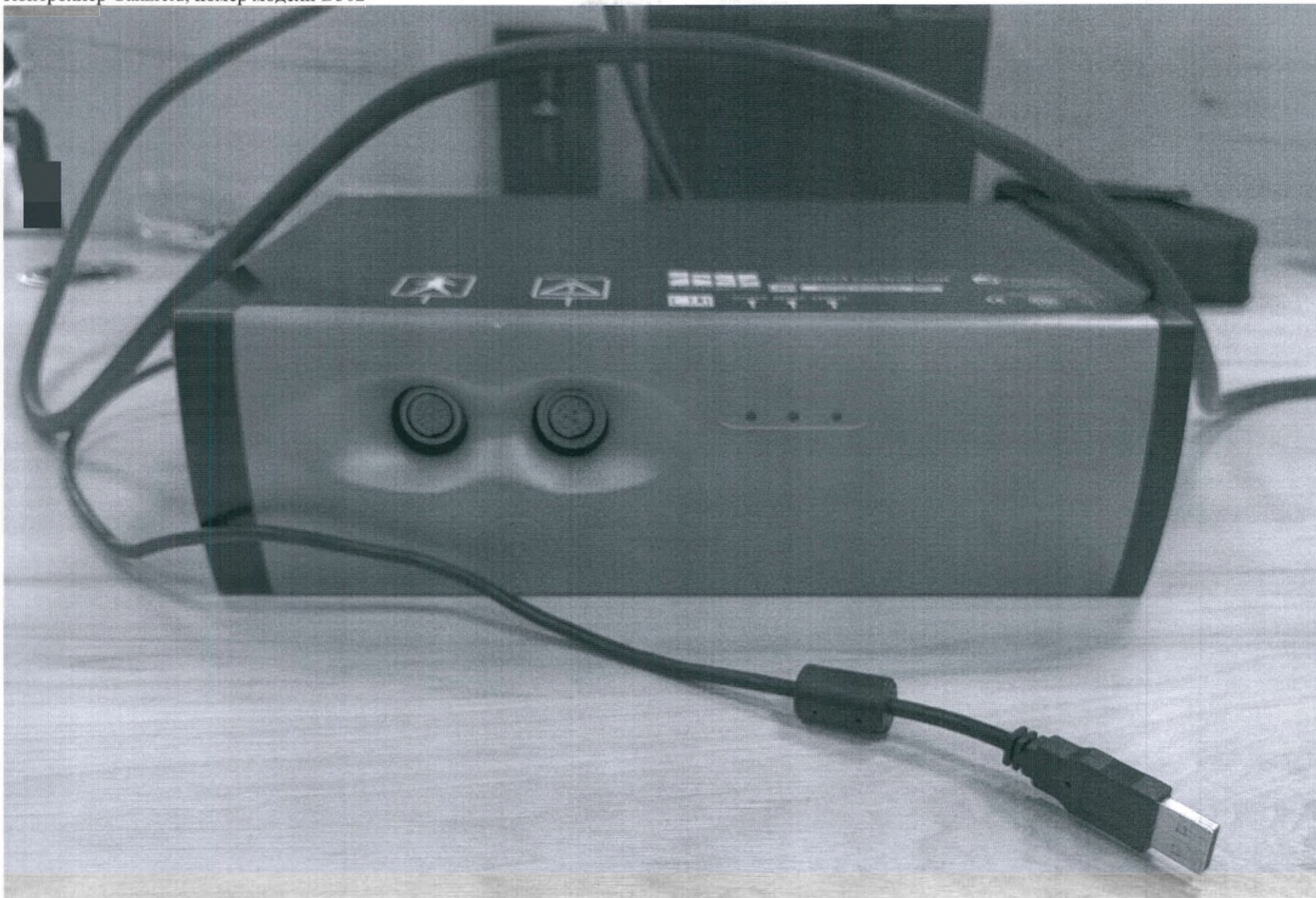








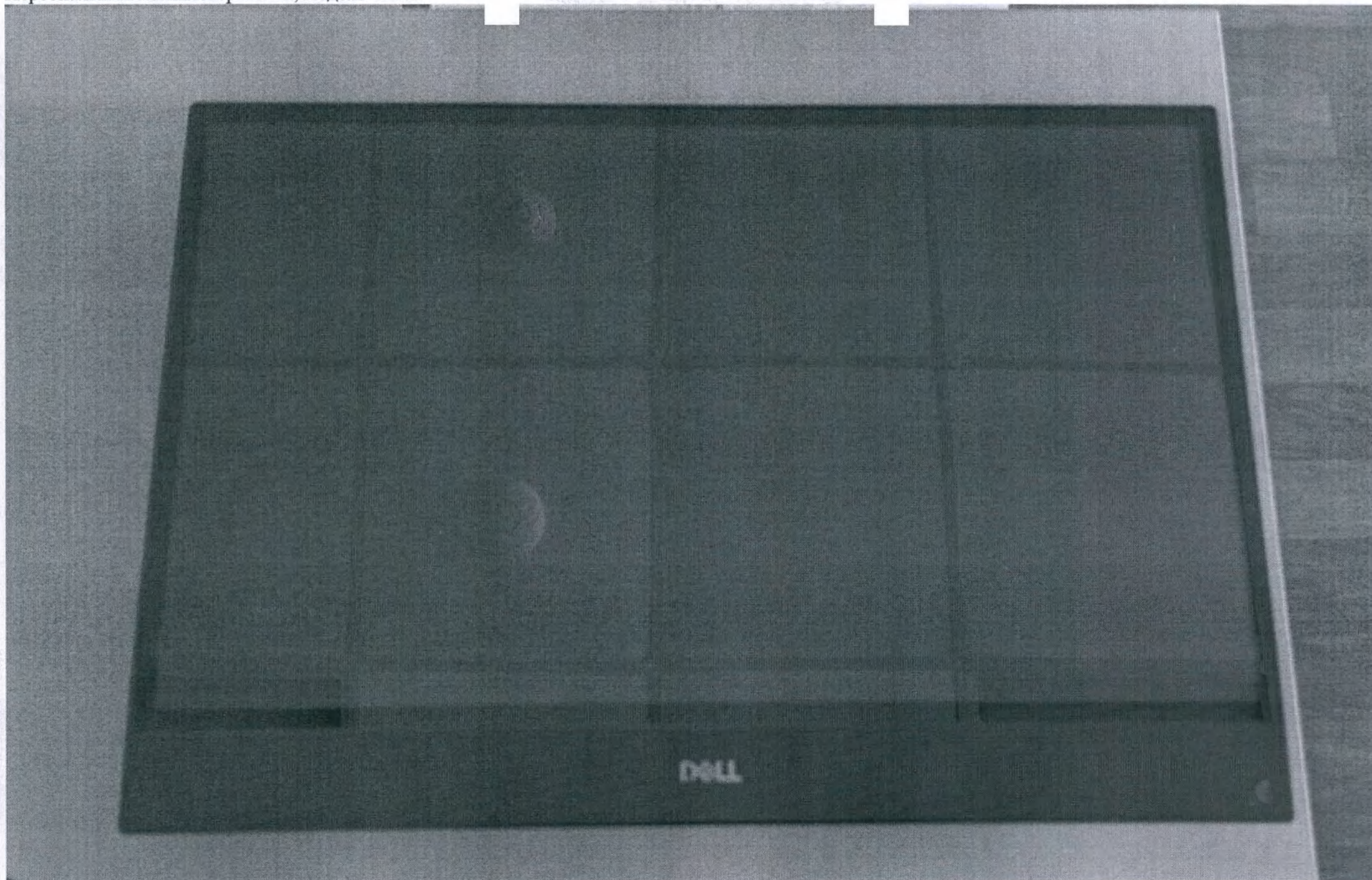


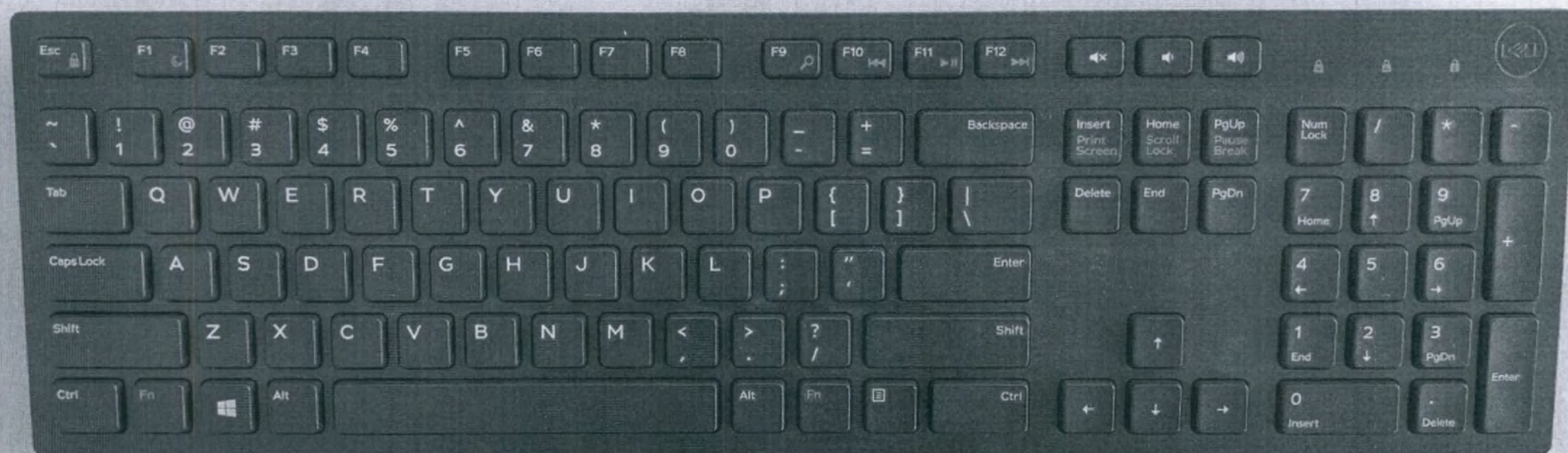

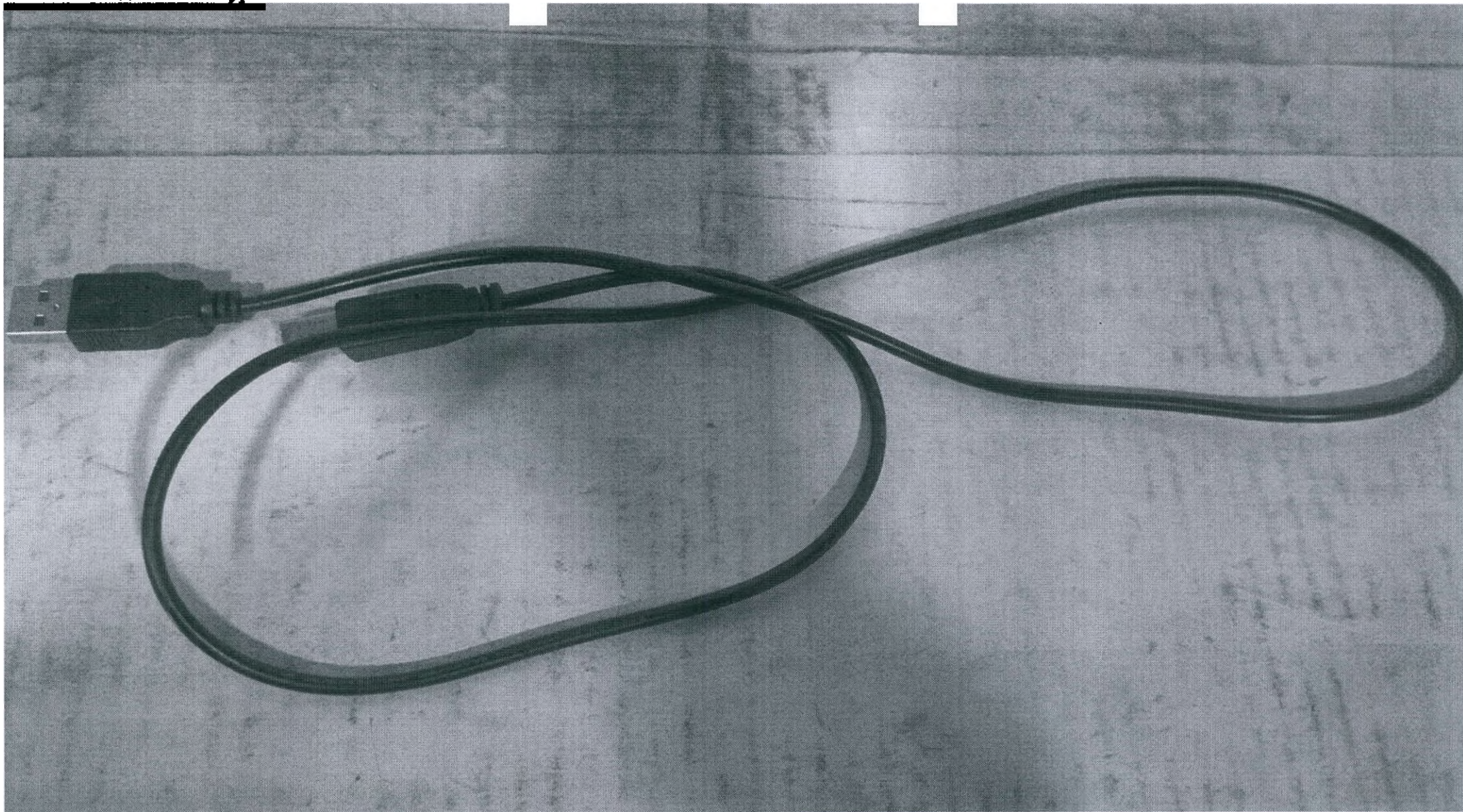


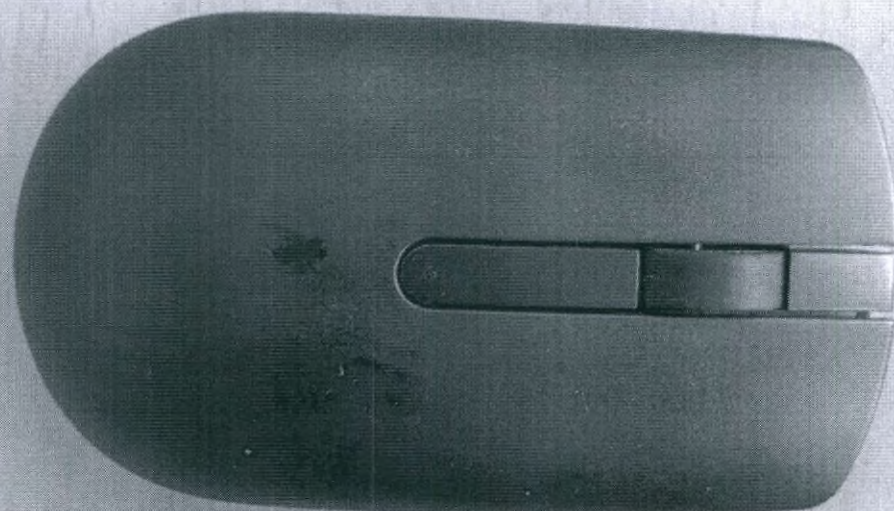


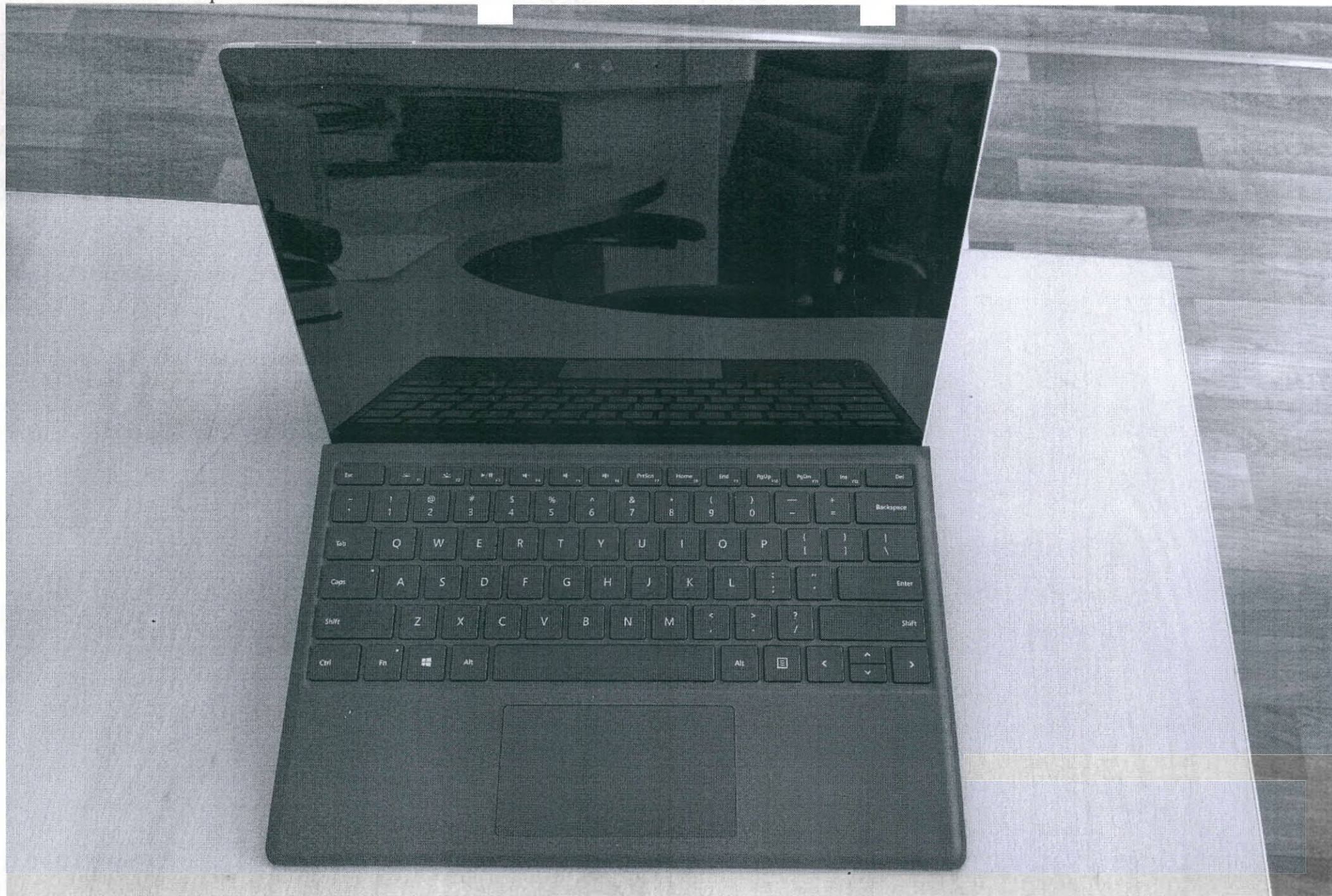


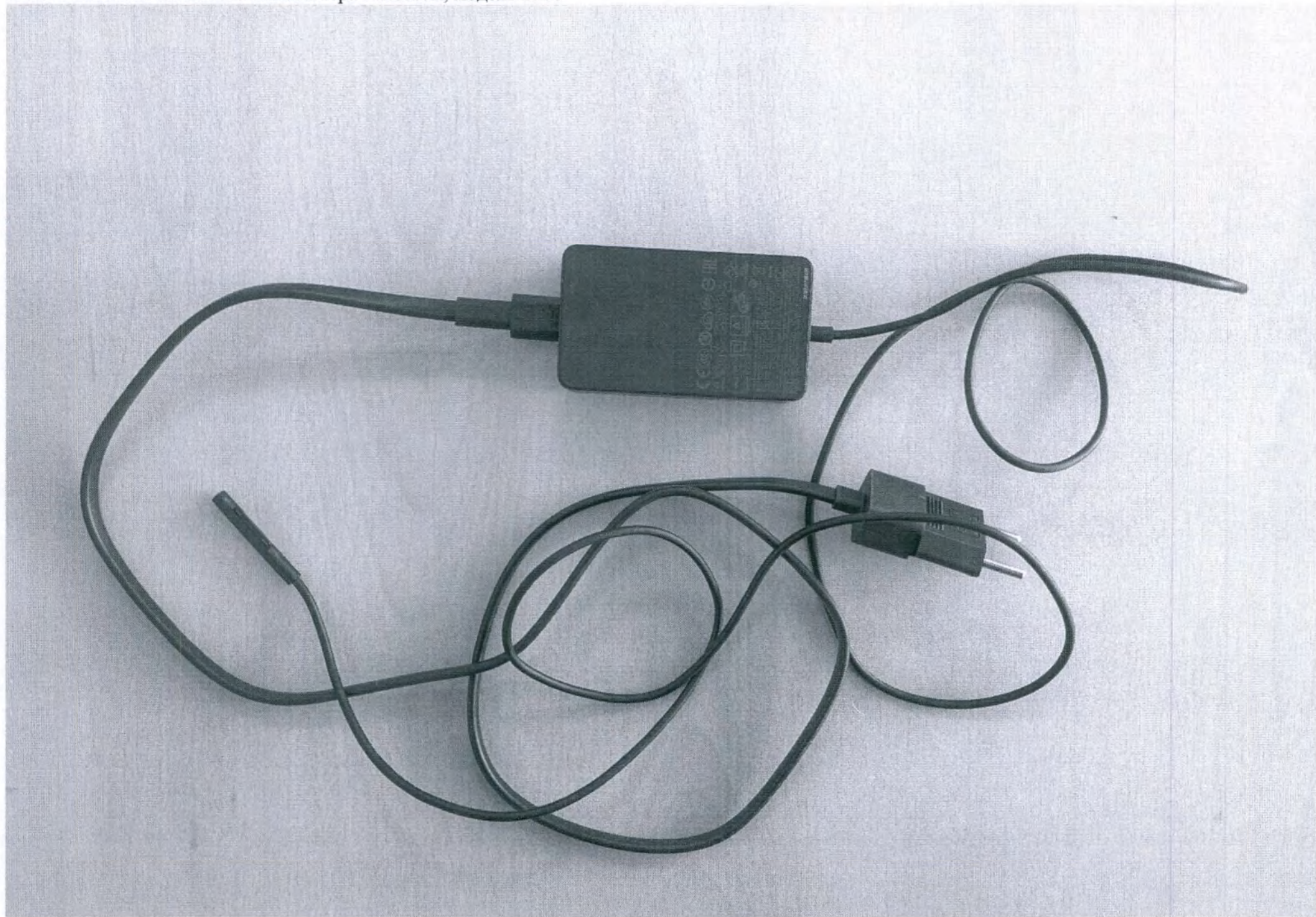




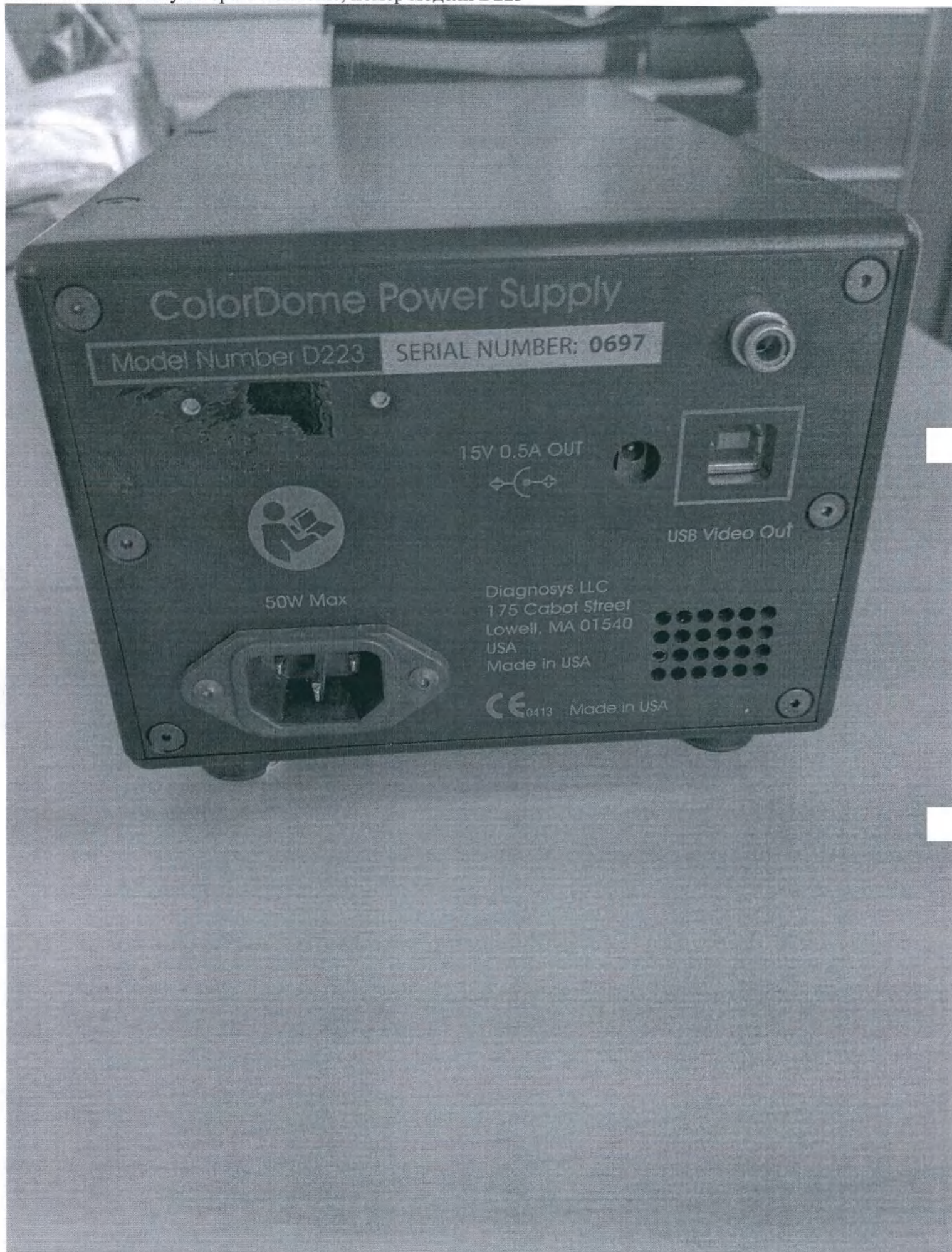

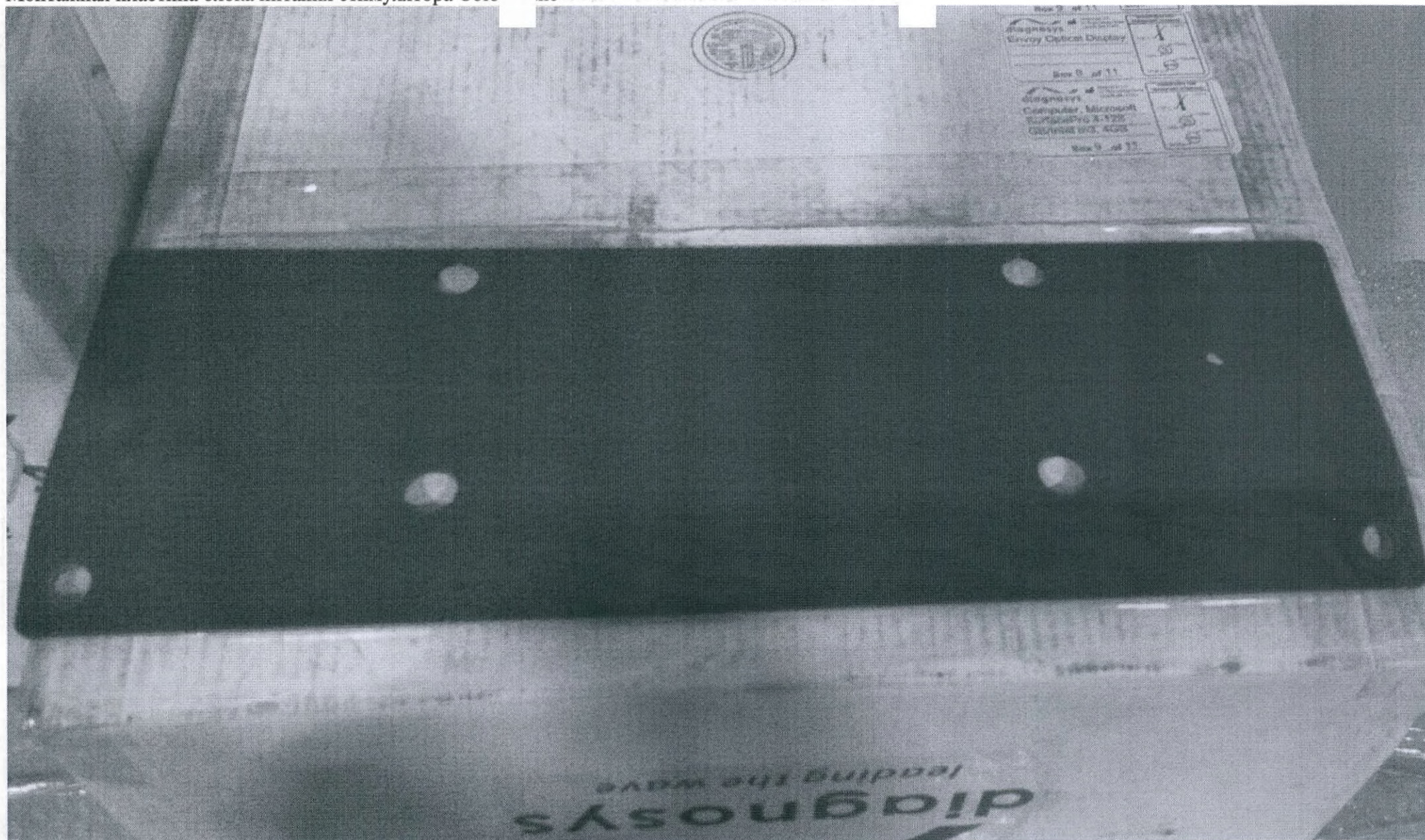


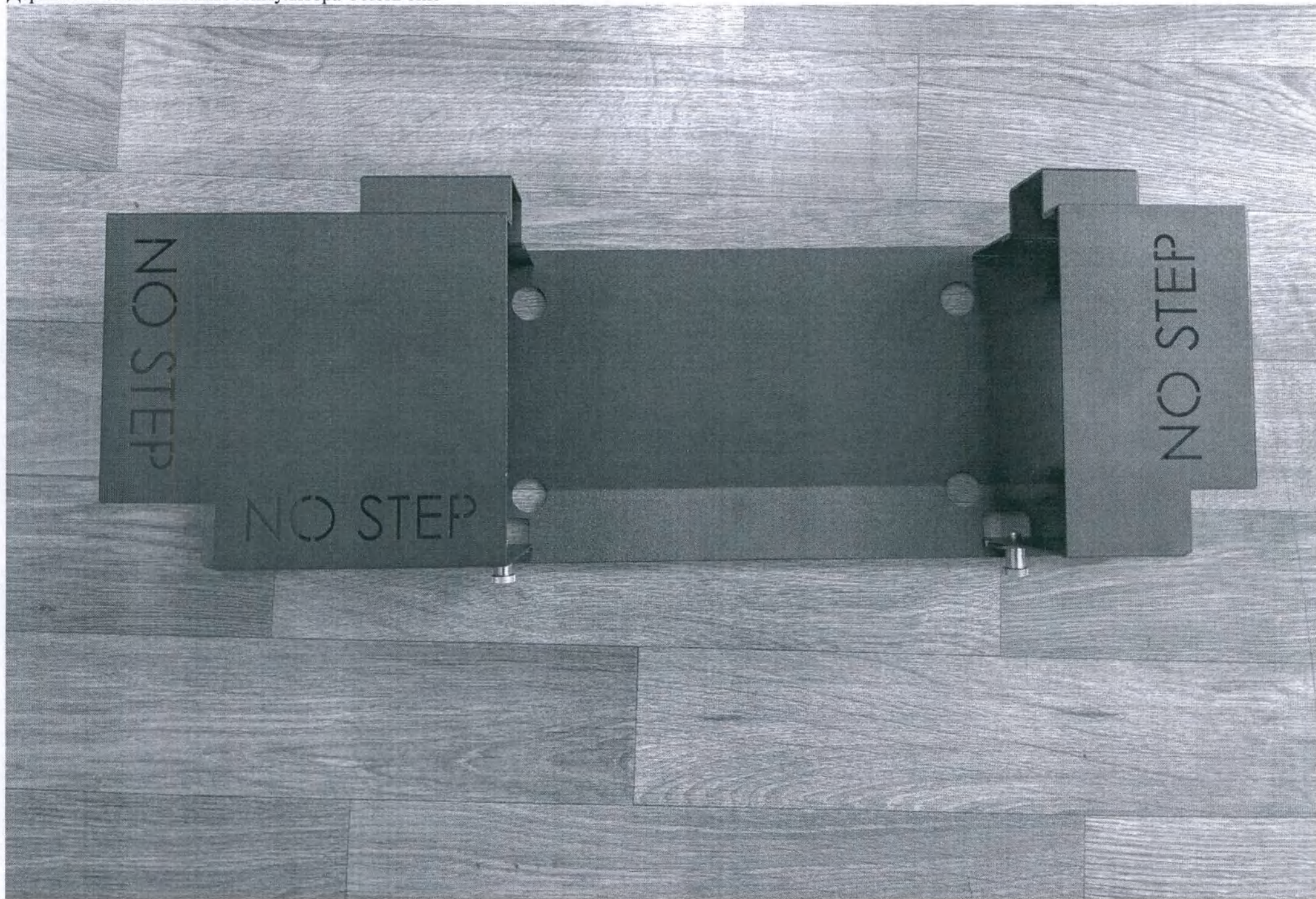





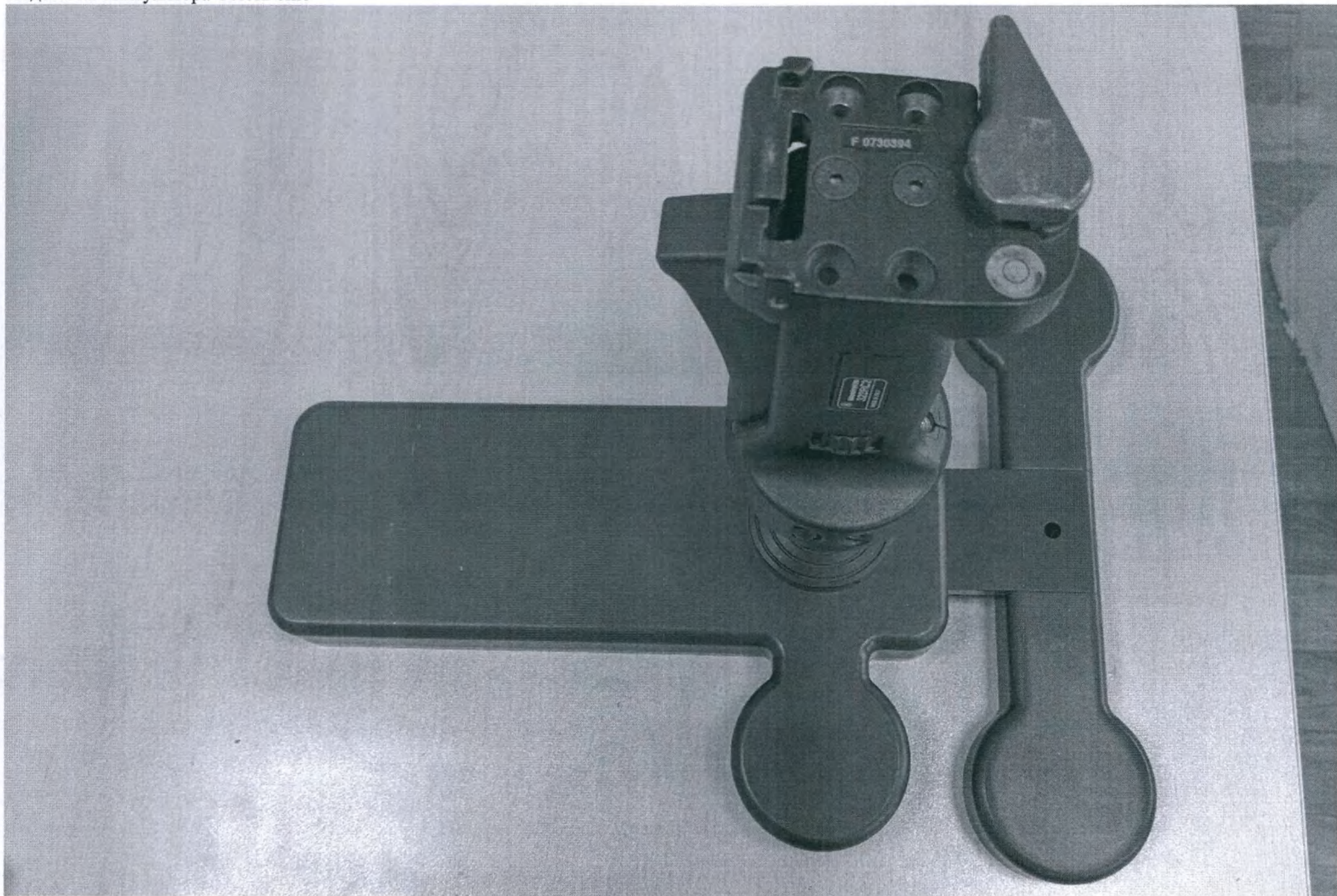




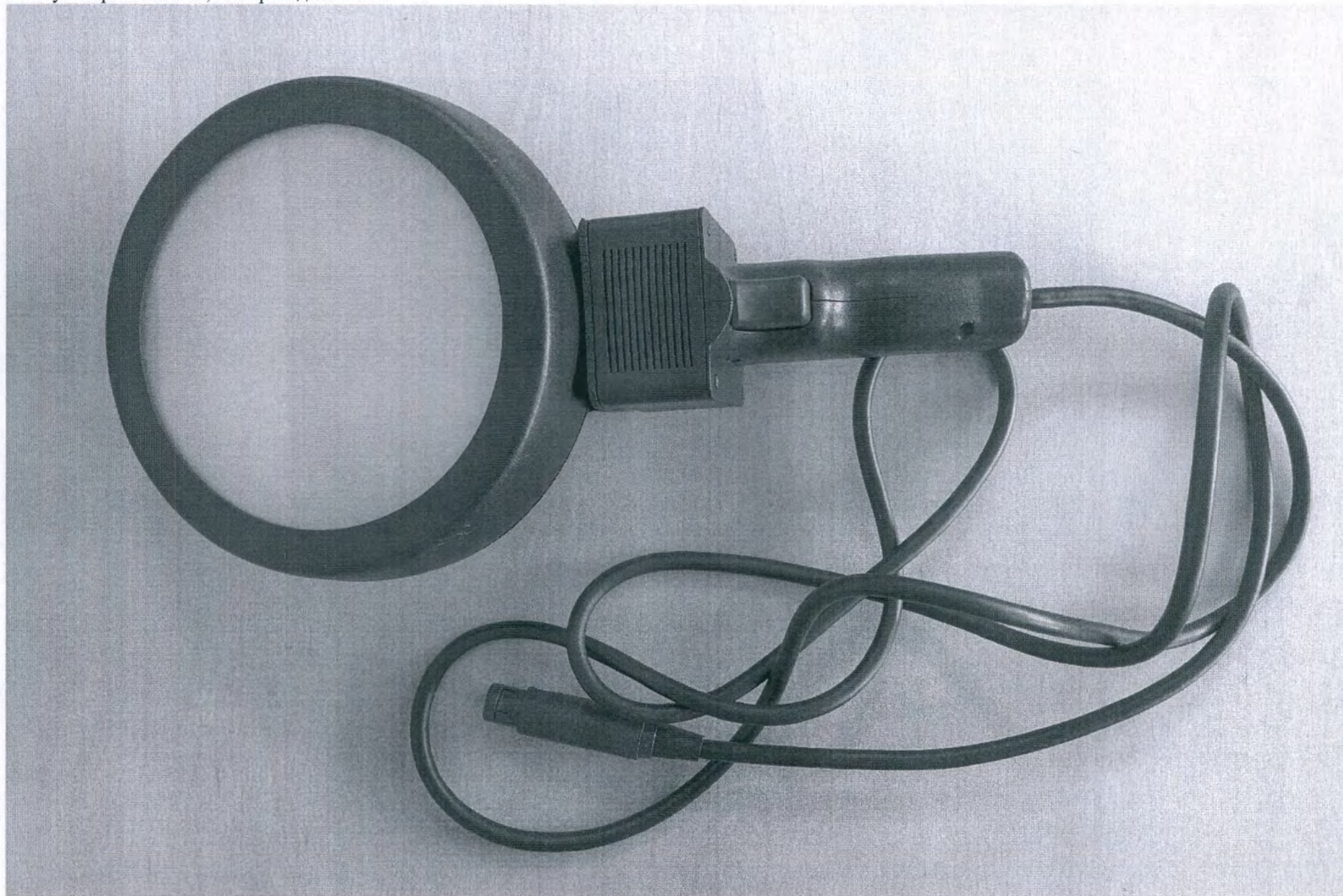

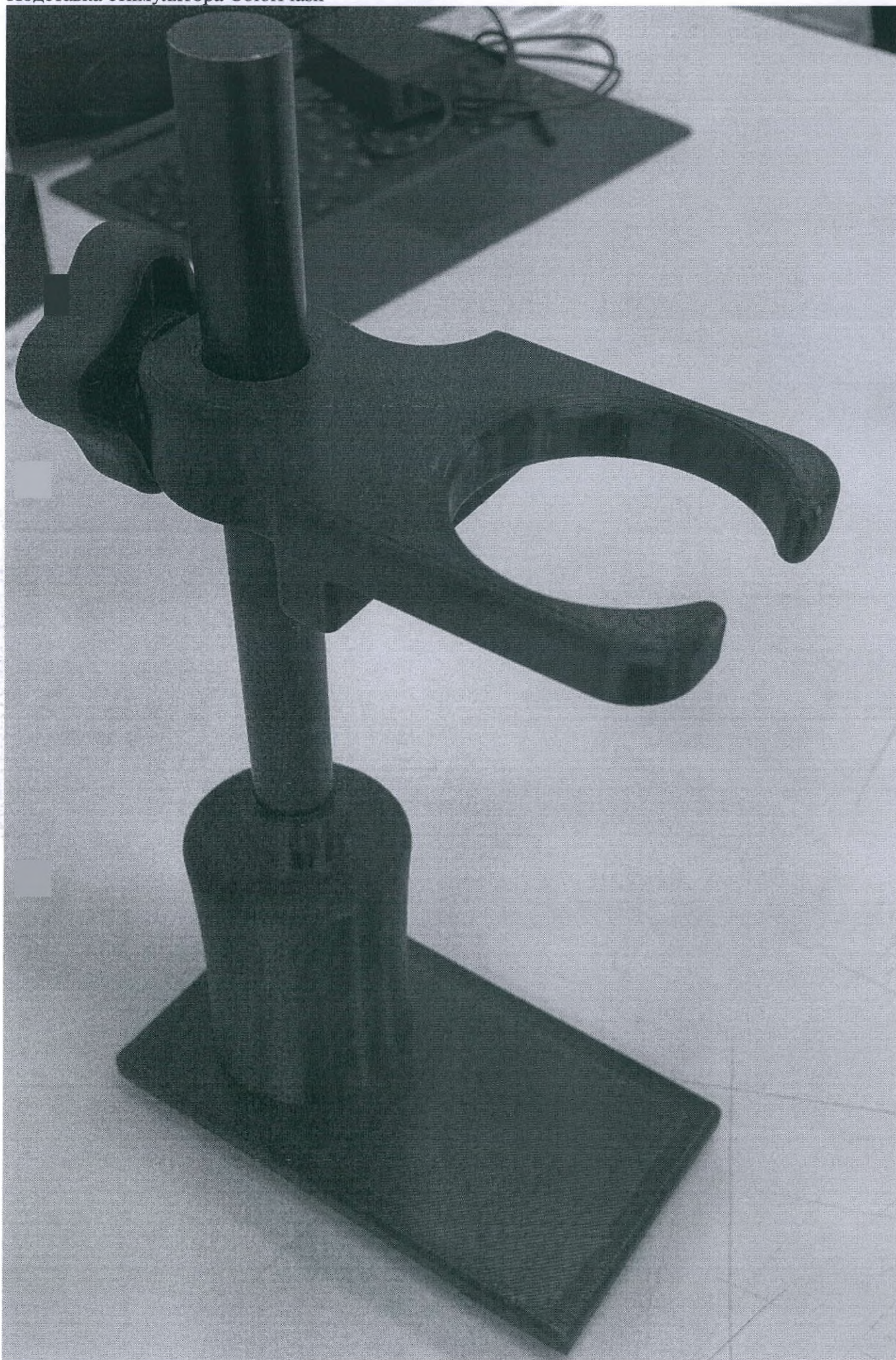


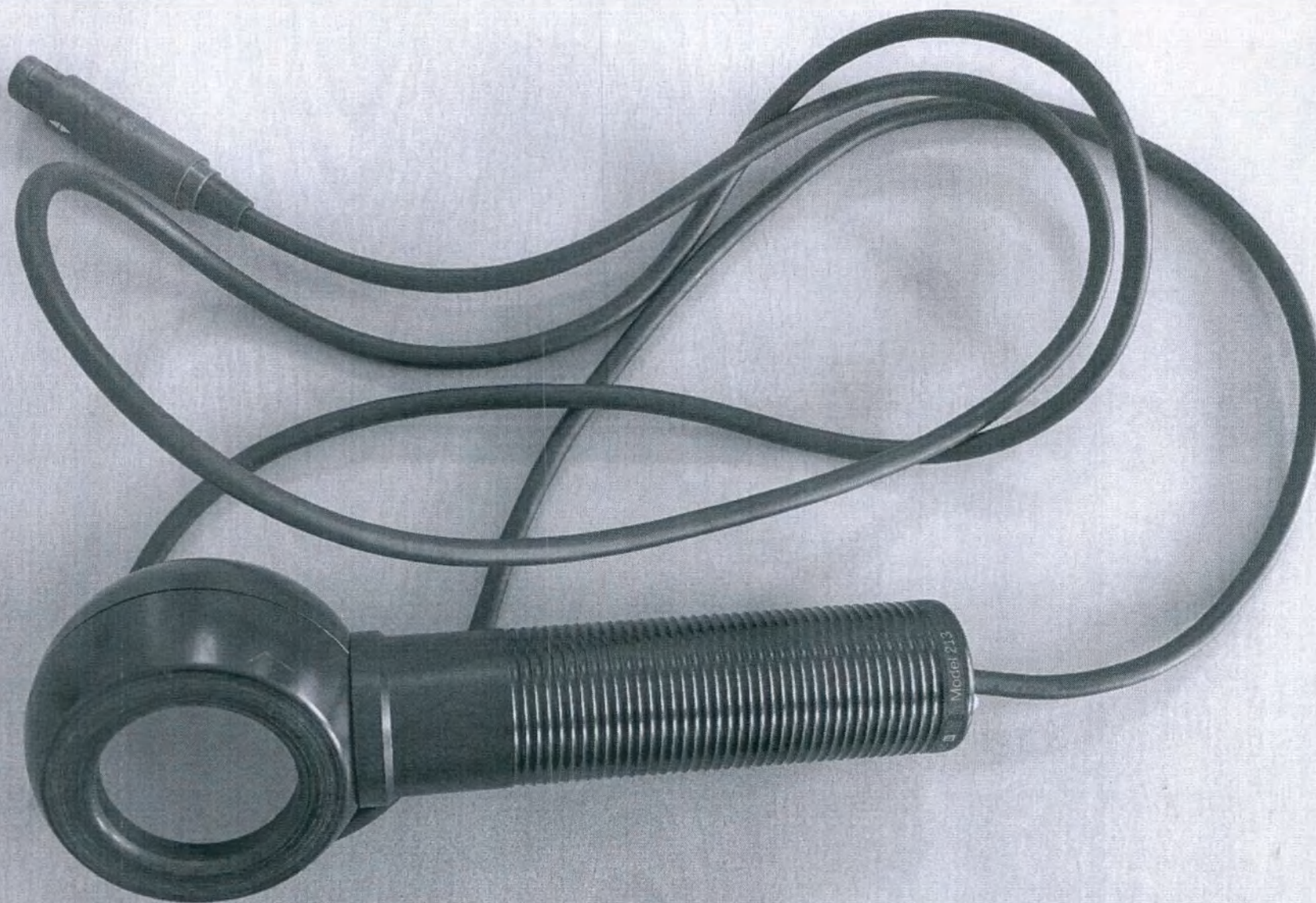


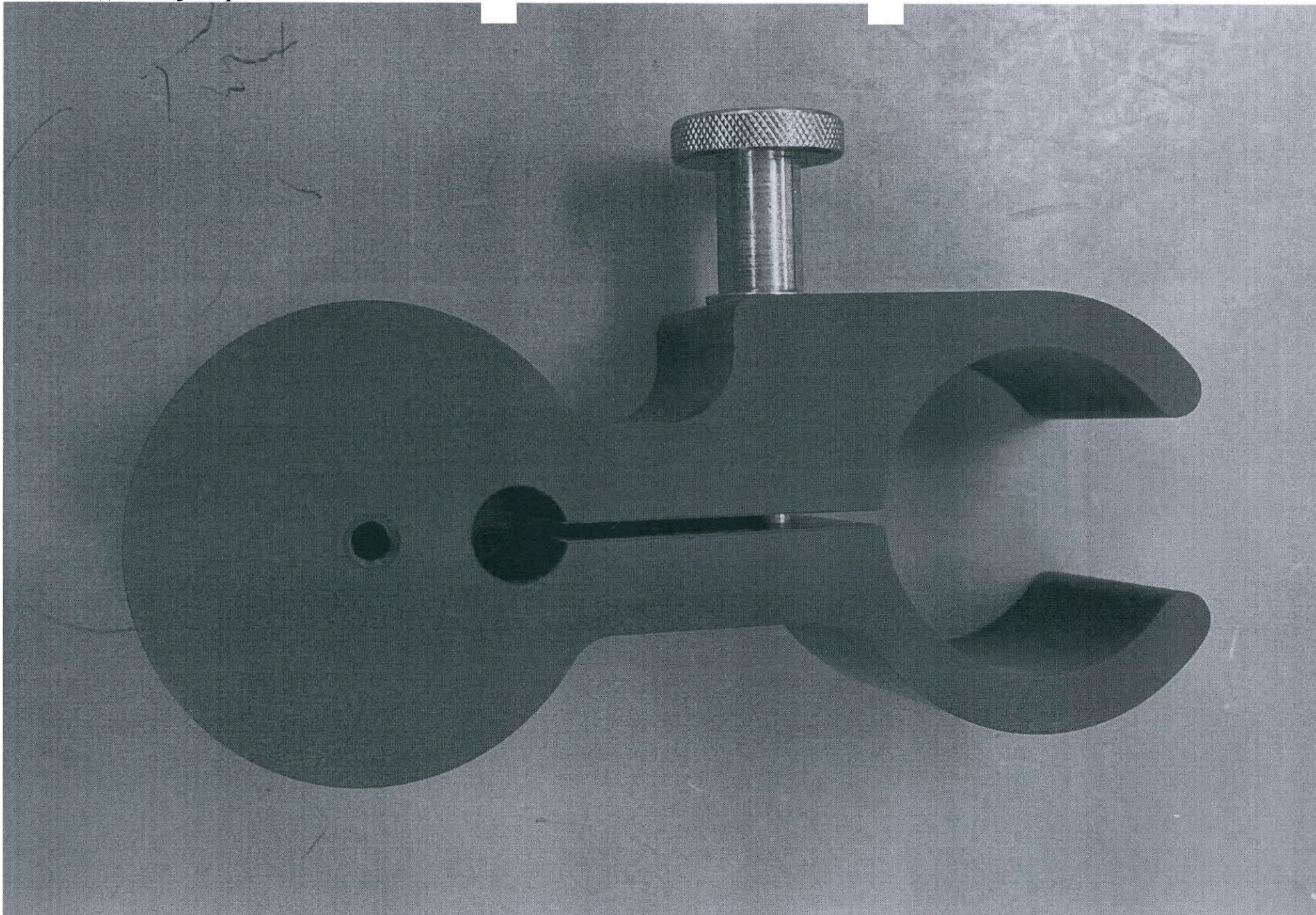







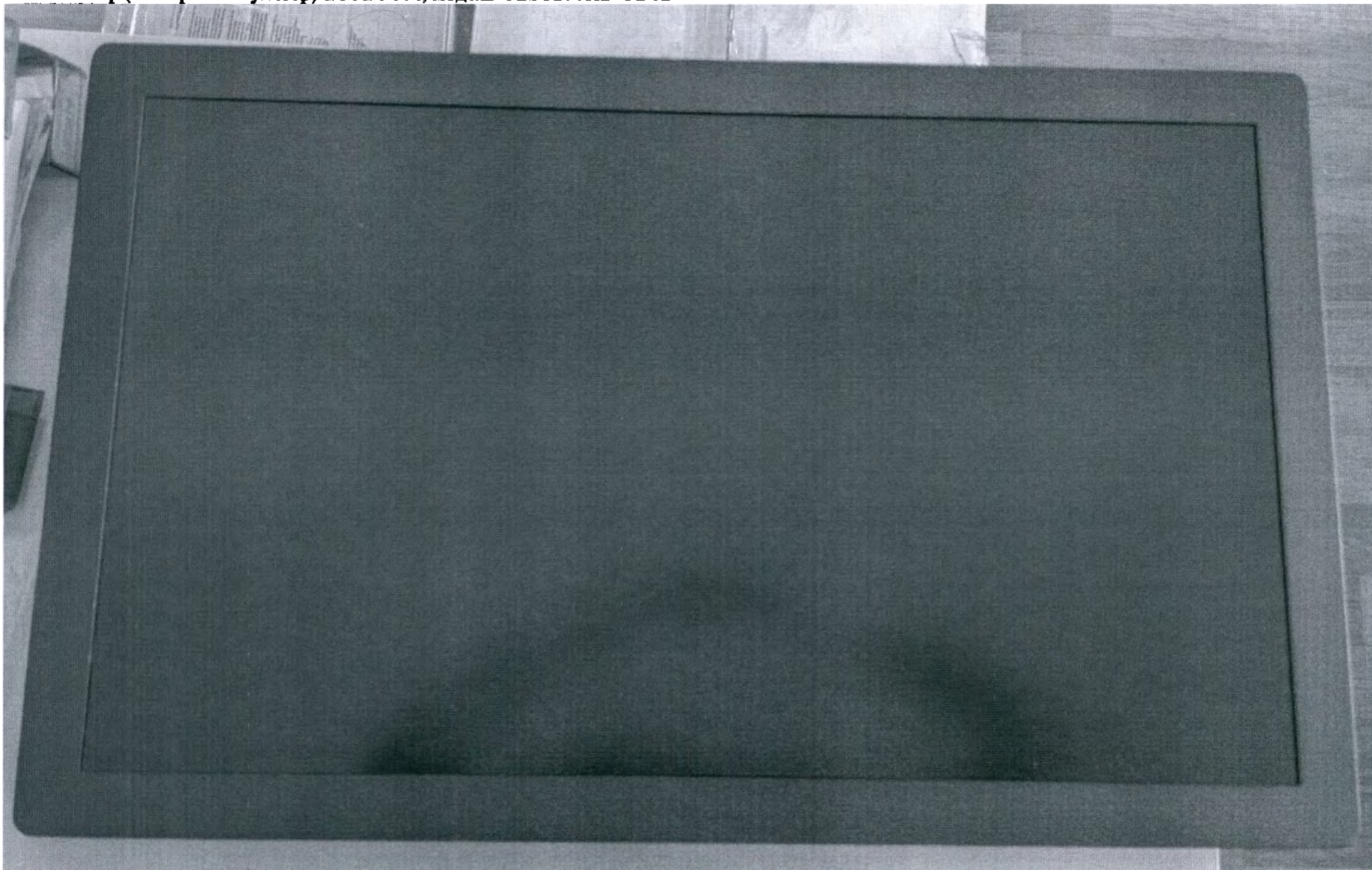


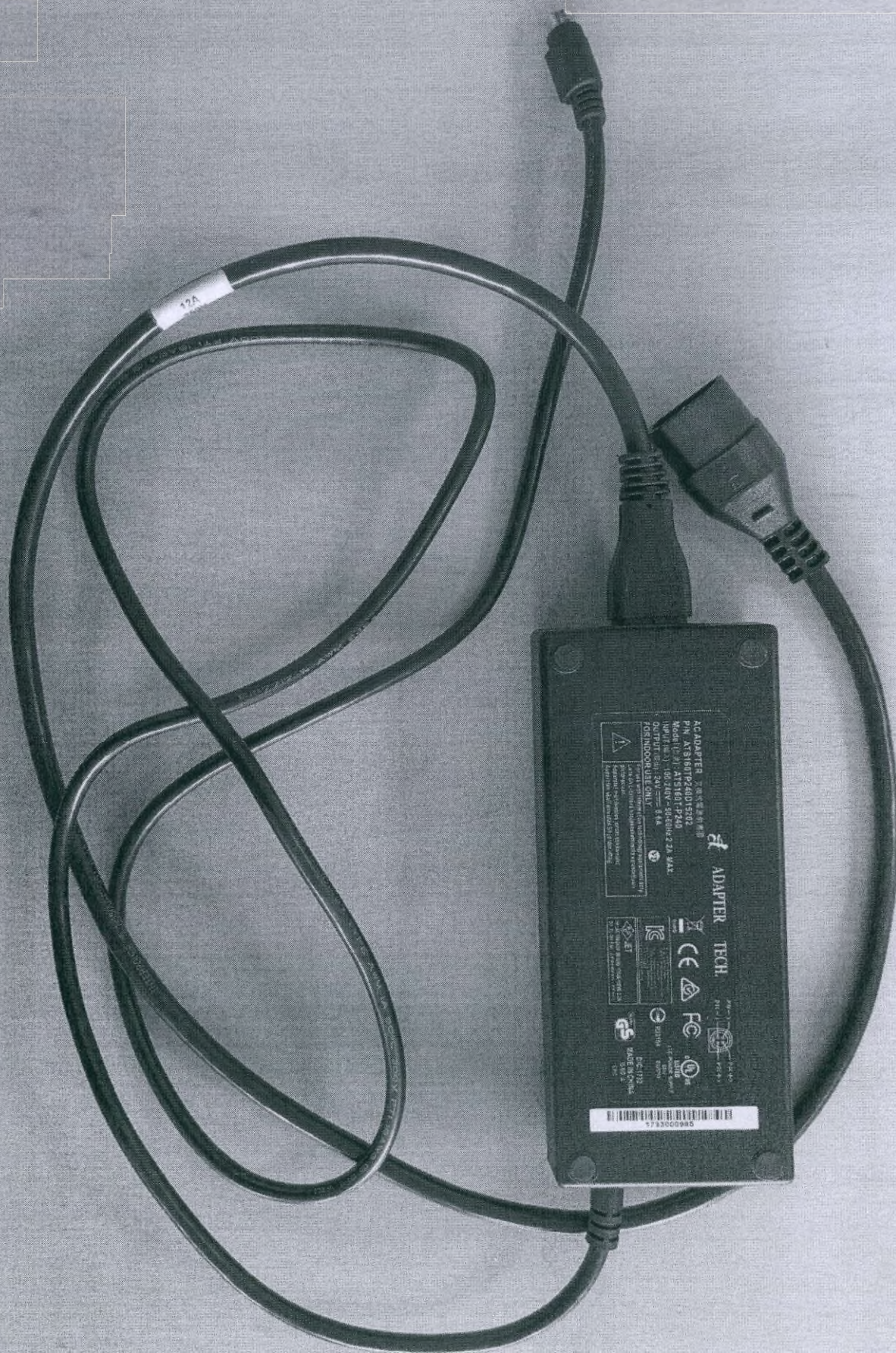

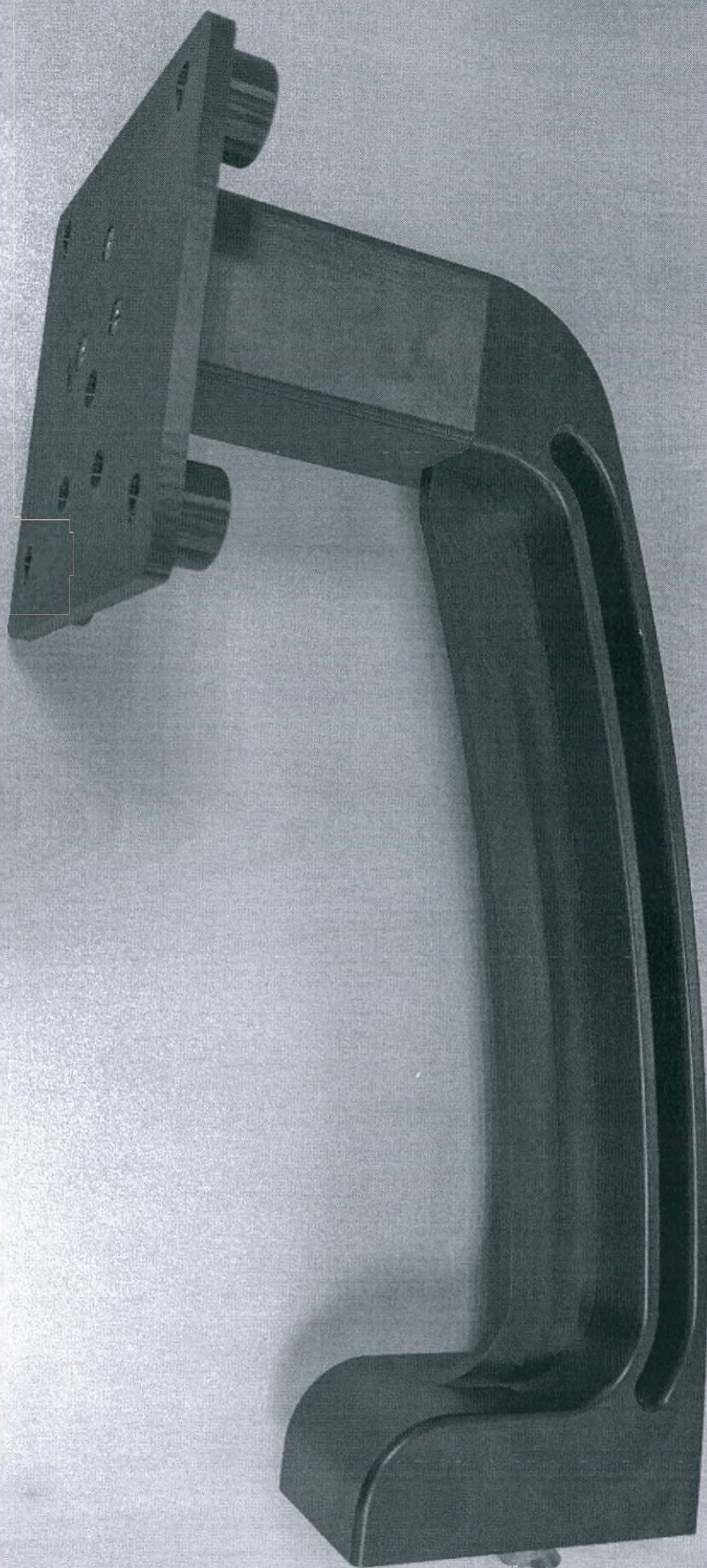


Подбородник в сборе

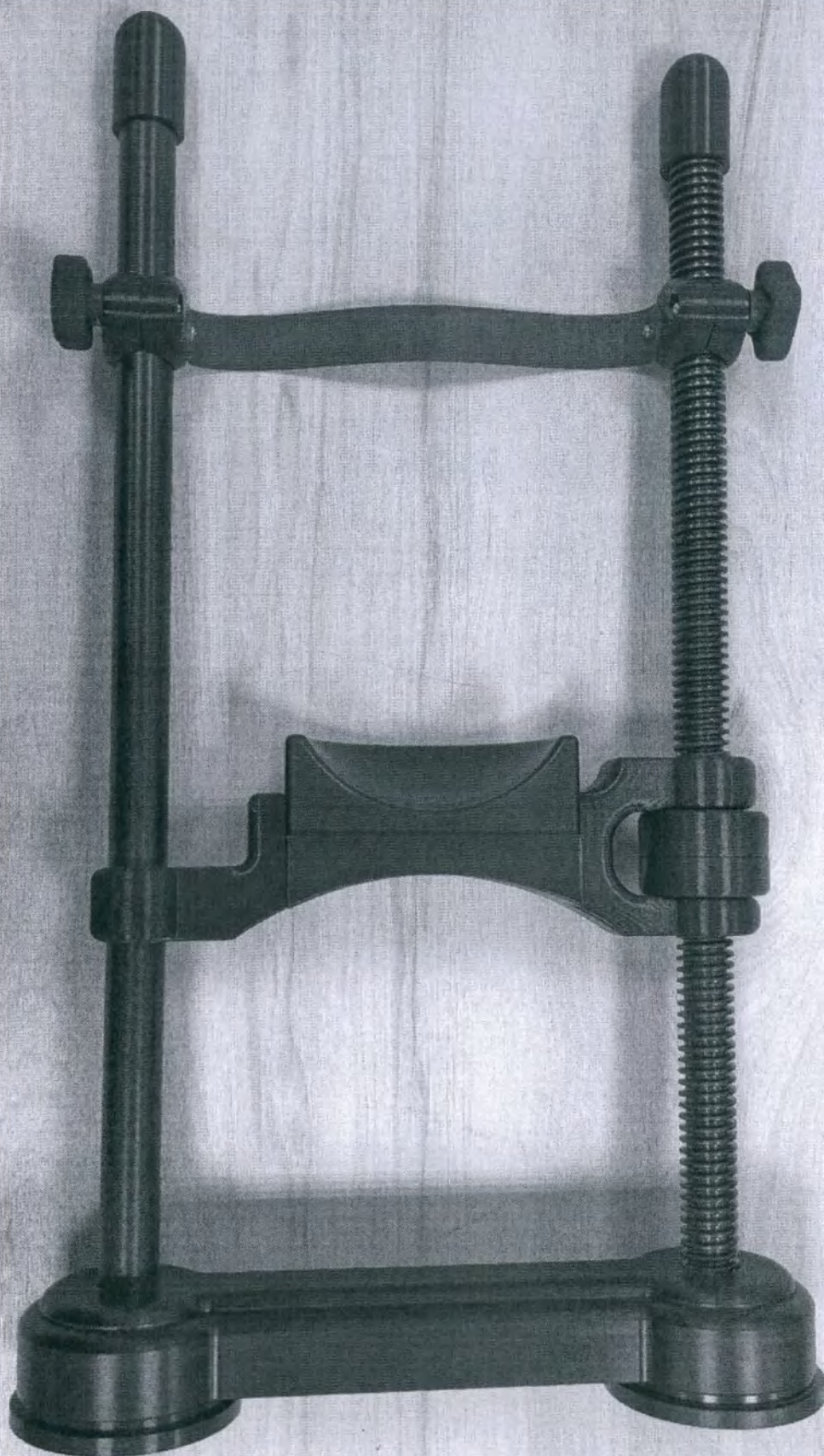

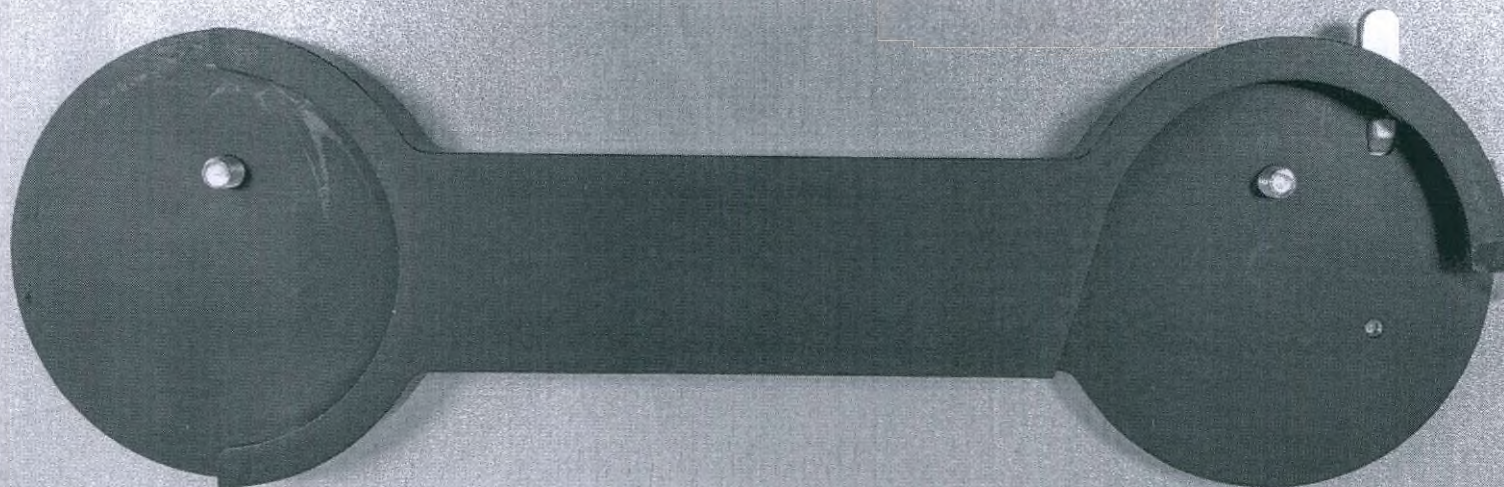

Кабели – не более 10 шт.

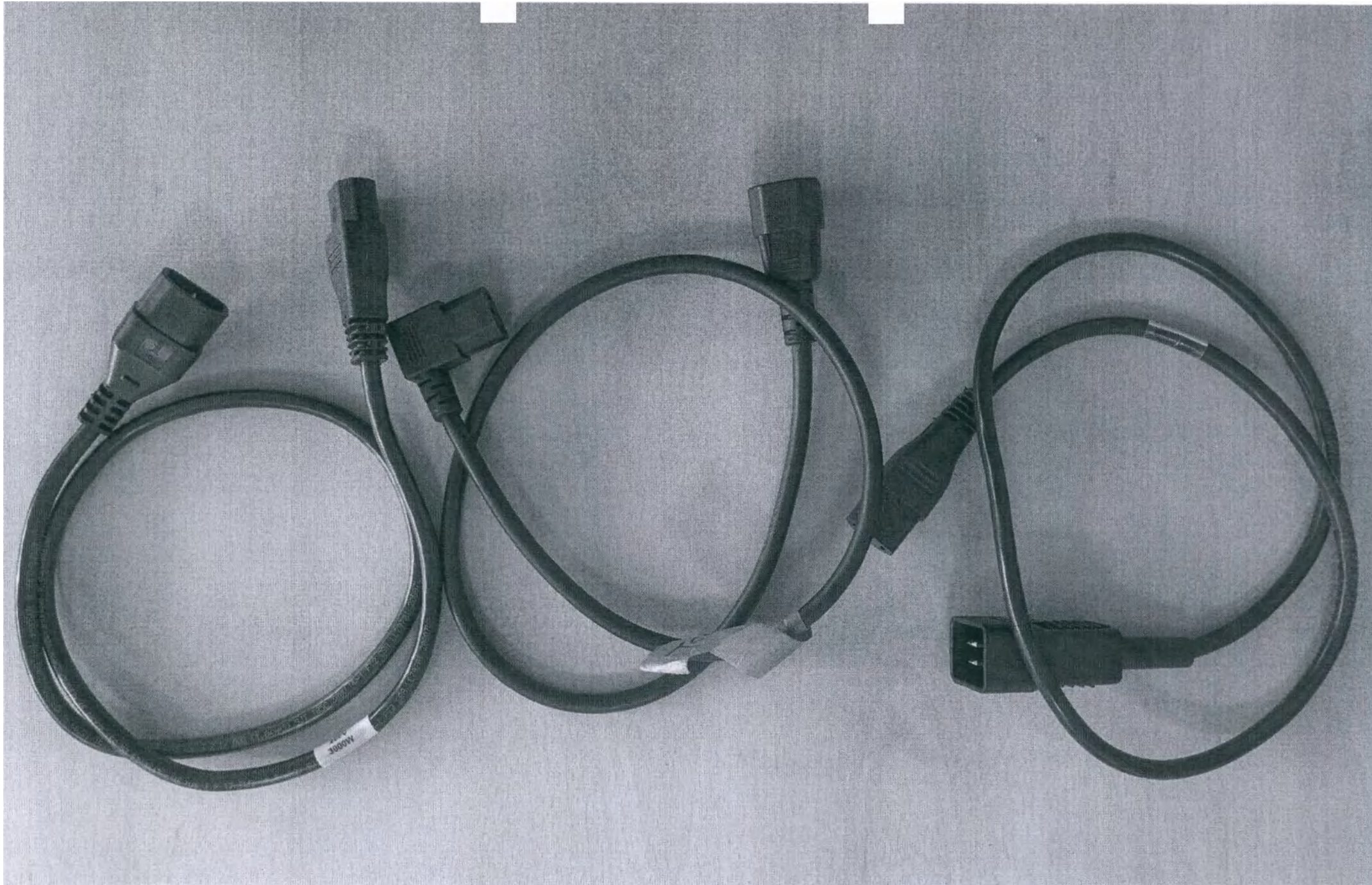



DIAGNOSYS DICOM: OVERVIEW


License Requirements

1. A DICOM license to enable the 3rd party software
2. A License code (not the same as 1) needs to be entered into the Espion System.

The Espion software requires the rzdcx.dll in either 32bit or 64 bit versions to be installed

Espion Setup

- In Configure System, "PROGRAM", "GLOBAL", is an option "DICOM License Code" which needs to be filled in. This is required before DICOM can be enabled
- In Configure System, under "PROGRAM", "GLOBAL" is an option called "DICOM Printing Enabled". This has three options: "Disabled", "DICOM+PRINT", "DICOM Only". These control which buttons on the print menu are enabled or disabled. "Disabled" disabled the DICOM button. "DICOM+PRINT" enables both the DICOM and Print buttons, and "DICOM Only" disables the Print button. For example, the following is "Disabled":

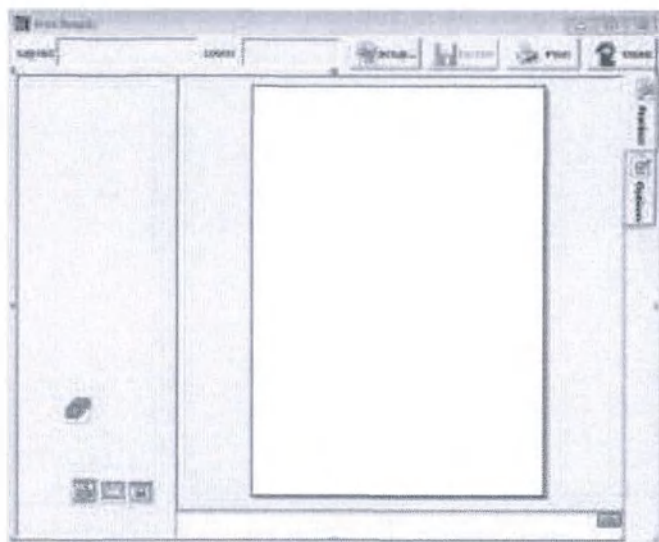

The DICOM button when pressed creates a DICOM compatible file with the extension ".DCM"

- In Configure System, under "PROGRAM", "GLOBAL" is an option called "DICOM Destination". This denotes the path (directory) where to save the DICOM files.

US: Diagnosys LLC USA • 175 Cabot Street, Suite 500 • Lowell, MA 01854 USA • Phone: 978.458.1600 • Fax: 978.458.1755 • Email: sales@diagnosysllc.com
UK: Diagnosys LTD UK • 5 Trust Court, Vision Park • Histon, Cambridge, CB24 9PW • Phone/Fax: +44 (0) 1223 520699 • Email: mail@diagnosysuk.co.uk

www.diagnosysllc.com

Total records: 115 [Name: TEST ID: tests: 3, last: 21-Dec-2021] [SYSDBA-ADMINISTRATOR-Hardware On : V11]

Method	Stimulus	Totals	Name	Sex	Age	Range	Comment	Electrode Type
		20	A Wave		0	150 yrs	SCOTOPIC RED	SKIN
		20	B Wave		0	150 yrs	SCOTOPIC RED	SKIN
		20	A Wave		0	150 yrs	SCOTOPIC BLUE	SKIN
		20	B Wave		0	150 yrs	SCOTOPIC BLUE	SKIN
		20	A Wave		0	150 yrs	SCOTOPIC WHITE	SKIN
		20	B Wave		0	150 yrs	SCOTOPIC WHITE	SKIN
		20	A Wave		0	150 yrs	SCOTOPIC WHT MAX	SKIN
		20	B Wave		0	150 yrs	SCOTOPIC WHT MAX	SKIN
		20	A Wave		0	150 yrs	PHOTOPIC WHITE	SKIN
		20	B Wave		0	150 yrs	PHOTOPIC WHITE	SKIN

Totals: Name Filter: Sex Filter:

New
 Modify
 Delete
 Transfer
 Copy
 First
 Last
 Prior
 Next
 Close


Total records: 3 [Name: TEST ID: tests: 3, last: 21-Dec-2021] [SYSDBA-ADMINISTRATOR-Hardware On : V11]

Method	Stimulus	Totals	Name	Sex	Age	Range	Comment	Electrode Type
		16	All Ages		0	40 yrs	A861-404C-B7BA-B739	DTL Fiber
		106	All Ages		40 yrs	20 yrs	A861-404C-B7BA-B739	DTL Fiber
		130	All Ages		60 yrs	150 yrs	A861-404C-B7BA-B739	DTL Fiber

Totals: Name Filter: Sex Filter:

New
 Modify
 Delete
 Transfer
 Copy
 First
 Last
 Prior
 Next
 Close

Модуль DiagnosysFST:
- Кнопочный пульт, номер модели D196

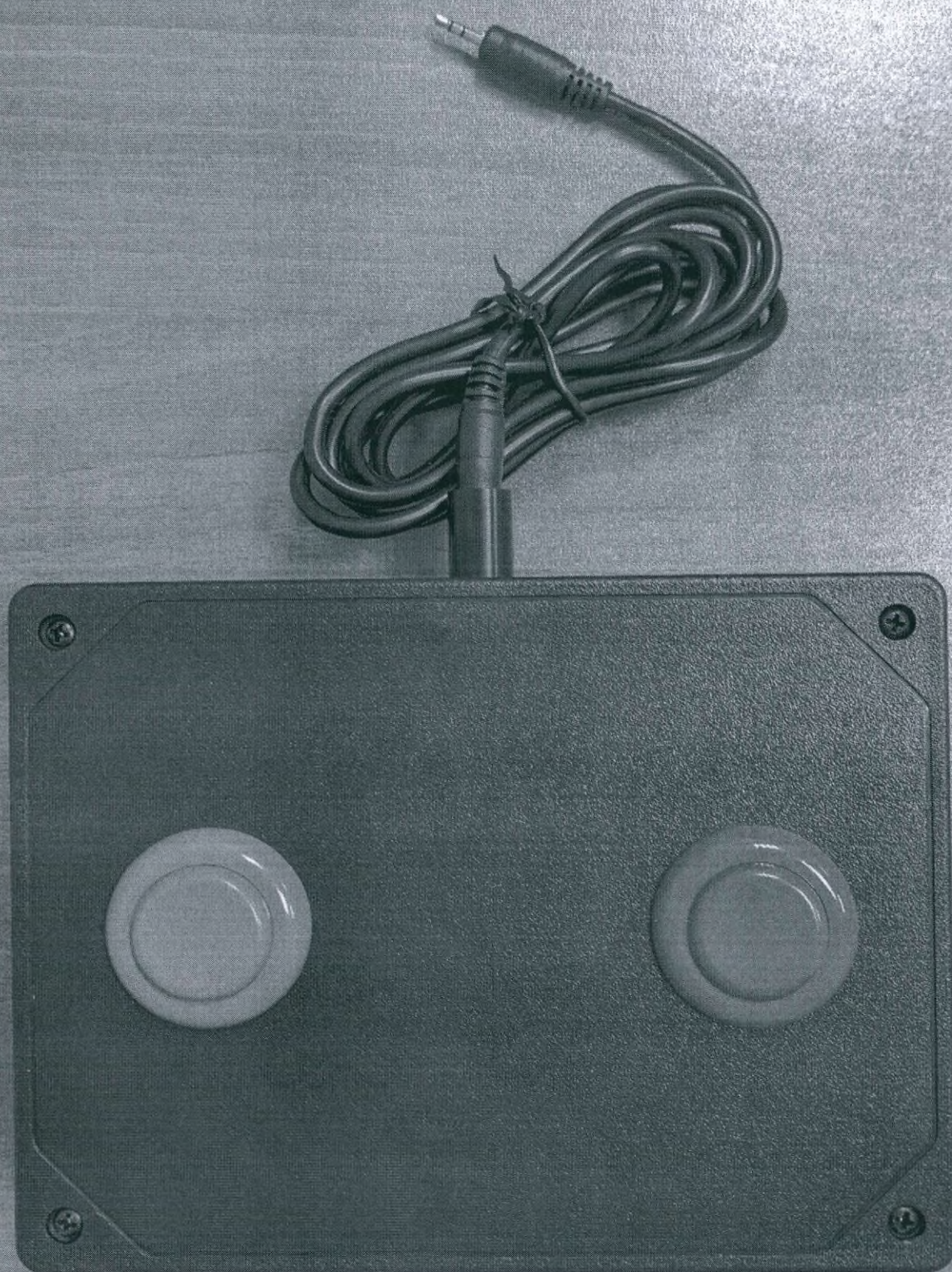

DiagnosysFST®

Module: License Key



DiagnosysFST® Module License

License Key: xxx-xxxxxxx

License holder:

Name: _____

Address: _____

Date: _____

Fair Use License & Limited Warranty

1. We grant you one license to install and use this software on a single computer at a time. If you do not agree to the following terms of this license, please uninstall and remove all copies and return the product to the place that you purchased it from within 30 days of your purchase for a full refund.
2. You may install and use the FST module on another computer, but the FST module may not be in use on more than one computer at a time unless you purchase additional licenses. You may make back-up copies of the software for archival purposes. You may permanently transfer your license to use the software to another party who will be bound by this agreement, provided you do not retain any copies of the software.
3. The software is protected by the copyright laws of the U.S. and other countries, and we retain all intellectual property rights in the software. You may not separately publish, sell, market, distribute, lend, lease, rent, or sublicense the software code. However, this license is not to be construed as prohibiting or limiting any fair use sanctioned by copyright law, such as permitted library and classroom usage.
4. We warrant that the software will provide the features and functions generally described in the product specification on our web site when you purchased it and in the product documentation. Media on which the Software is furnished, if any, will be free from defects in materials and workmanship.
5. We have taken all reasonable steps to keep the software free of viruses, spy ware, "back door" entrances, or any other harmful code. We will not track or collect any information about you, your data, or your use of the software except as you specifically authorize. The software will not download or install patches, upgrades, or any third party software without getting your permission. We will not intentionally deprive you of your access to your data.
6. We do not warrant that the software or your ability to use it will be uninterrupted or error-free. To the extent permitted by applicable law, we disclaim any implied warranty of merchantability or fitness for a particular purpose.
7. Your exclusive remedy under the above limited warranty shall be, at our option, either a full refund of the purchase price or correction of the defective software or media. To the fullest extent permitted by applicable law, we disclaim all liability for indirect or consequential damages that arise under this license agreement. Nothing in this agreement limits our liability to you in the event of death or personal injury resulting from gross negligence, fraud, or knowing misrepresentation on our part.
8. If any part of this agreement is found to be invalid or unenforceable, the remaining terms will stay in effect. This agreement does not prejudice the statutory rights of any party dealing as a consumer.
9. This agreement is governed by the laws, including Article 2 of the Uniform Commercial Code, of Massachusetts.
10. This agreement does not supersede any express warranties we made to you. Any modification to this agreement must be agreed to in writing by both parties.

www.diagnosysllc.com

US: Diagnosys LLC; 55 Technology Drive, Suite 100, Lowell, MA 01851; 978-456-1600; sales@diagnosysllc.com
EU: Diagnosys Vision Ltd; Office 117, DDC Building, Belheary Road, Swords, Dublin, K67 E5A0, Ireland; +44 (0) 1223 520699; mail@diagnosysvision.com
UK: Diagnosys UK Ltd; 5 Trust Court, Chivers Way, Vision Park Histon, Cambridge, CB24 9PW, UK; +44 (0) 1223 520699; mail@diagnosysuk.co.uk


The information listed herein is accurate as of the date of printing, however may change at any time without notice. The contents may differ from the current status of approval of the product in your country. Please contact your local Diagnosys representative for more information. © 2003 by Diagnosys, LLC. All rights reserved.
Doc: 12474 Rev 6 RCH 1665 Date: 15 Feb 2003

- Протоколы тестирования DiagnosysFST

FST Full-field Sensitivity Test- [User:NOT SELECTED, Role:SYSDBA]

Patient	Age	Test Time	Status	Eye [LE]	Mask	Buttons	3:34 14
102 201	121 Years						

Individual Eye Graphs Both Eye Graph Setup

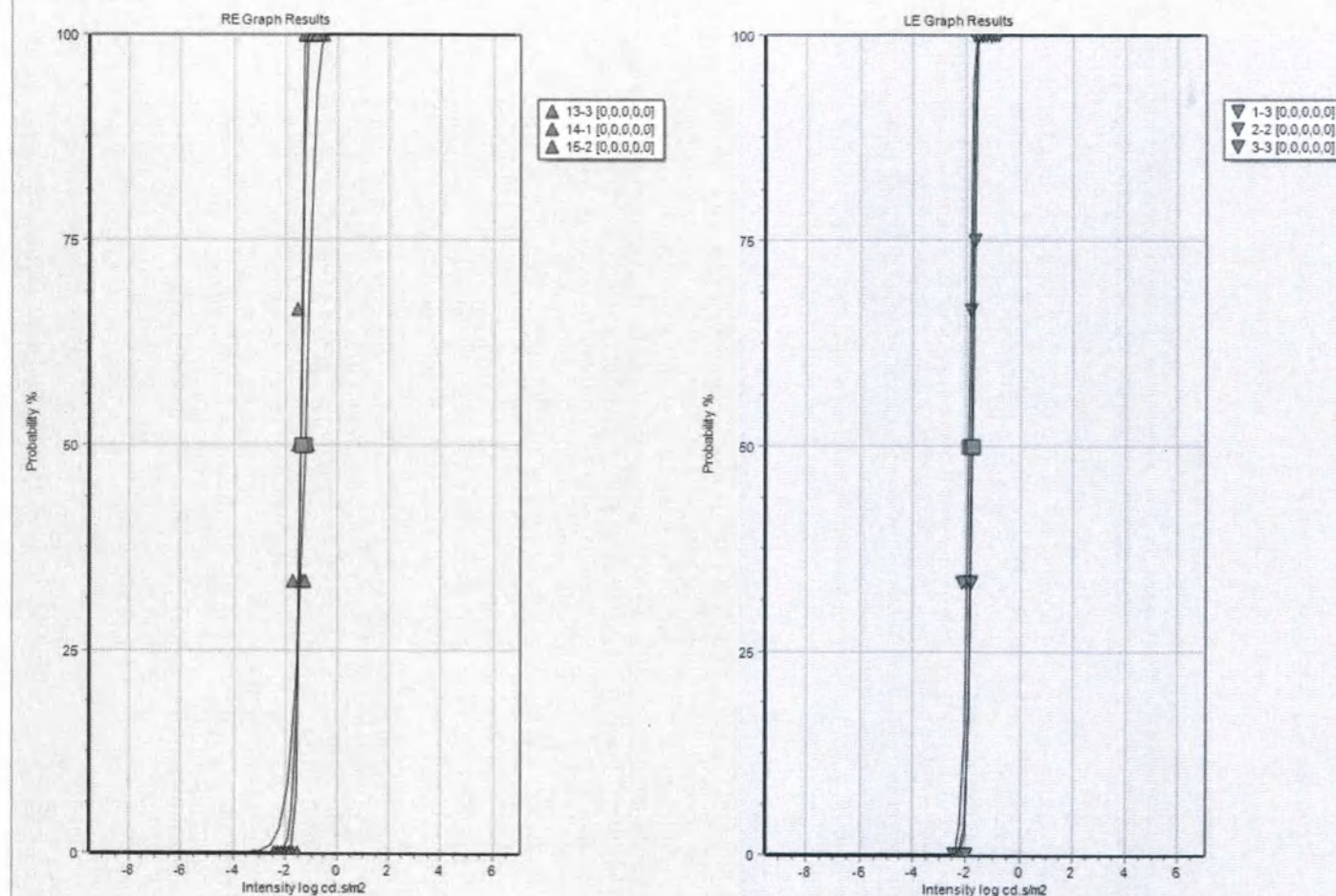

Results Run-time

All Test Results											
#	Color-FL	Eye	Combined Quality	Combined Threshold log cd.s/m2	EB%	EM%	TO	n	Period	n	Period
10	White	LE	3	-1.832	0	0	0	76	01:17		
22	White	RE	3	-1.357	11	0	0	90	01:29		

Individual Test Results											
#	Color-FL	Eye	Quality	Threshold log cd.s/m2	EB%	EM%	TO	n	Period	Time	Date
1	White	LE	3	-1.911	0	0	0	23	01:11	2:37:36 PM	10/28/2021
2	White	LE	2	-1.815	0	0	0	23	01:11	2:39:06 PM	10/28/2021
3	White	LE	3	-1.776	0	0	0	30	01:29	2:40:53 PM	10/28/2021
13	White	RE	3	-1.397	33	0	0	30	01:29	2:16:09 PM	10/28/2021
14	White	RE	1	-1.230	0	0	0	30	01:30	2:18:06 PM	10/28/2021
15	White	RE	2	-1.398	0	0	0	30	01:29	2:20:48 PM	10/28/2021

Display Options

Colors	Results	Raw Data	Curve Fit	Double Timeout	0dB->cd(s)/m2	cd(s)/m2 Units
White	Individual	☒	☒	☐	1	☒

Run

Practice

Stop

Pause

Test Eye

Timer

Locate

Align

Notes

Print

Export


Load

Quality

Exit

11°C Mostly cloudy

3:34 PM 1/25/2022

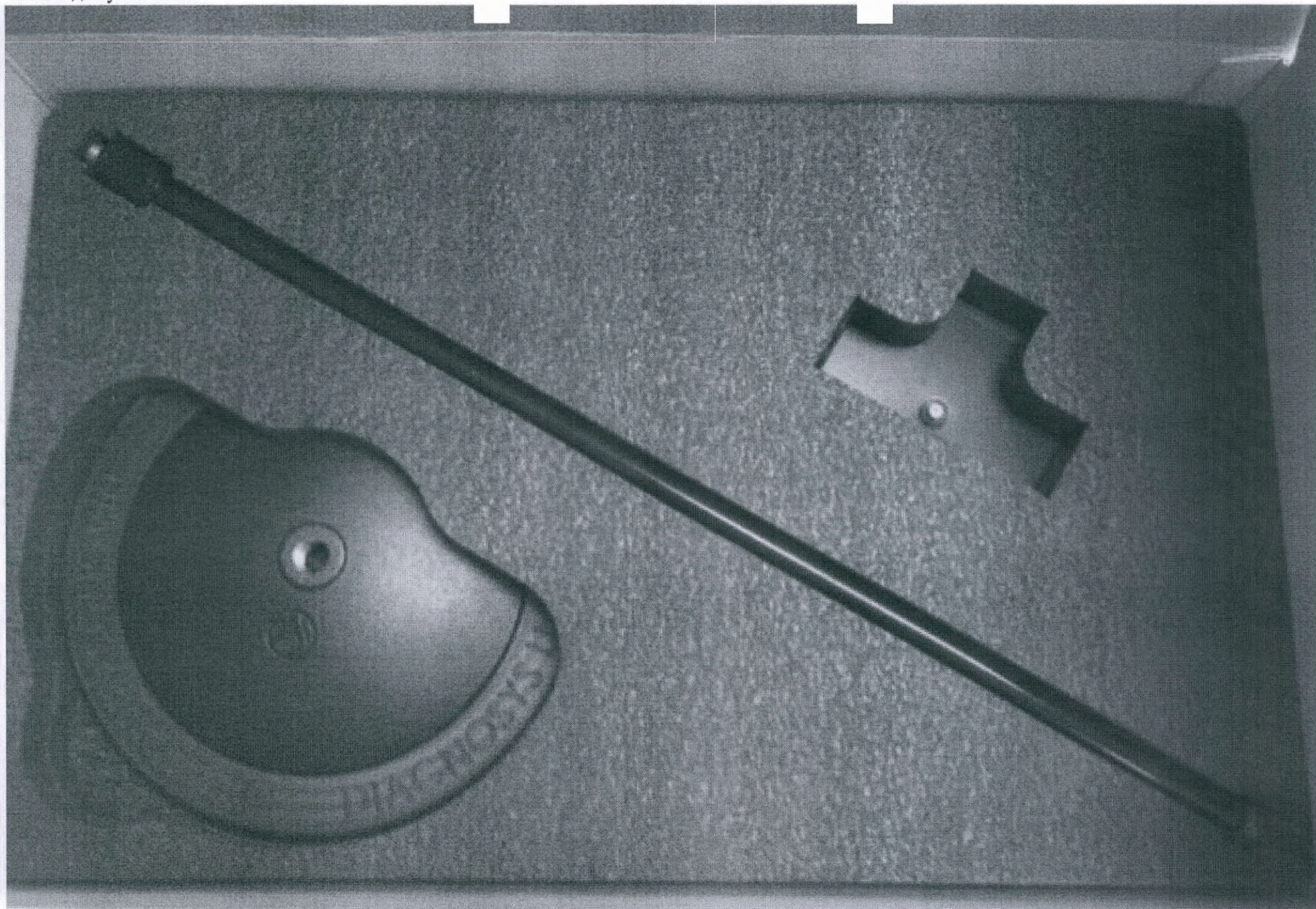

Модуль адаптации к темноте:

- Кнопка ответа пациента



Dark Adaptometry Module: License Key



Dark Adaptometry Module License

License Key: XXXX-XXXXXX

License holder:

Name: _____

Address: _____

Date: _____

Fair Use License & Limited Warranty

1. We grant you one license to install and use this software on a single computer at a time. If you do not agree to the following terms of this license, please uninstall and remove all copies and return the product to the place that you purchased it from within 30 days of your purchase for a full refund.
2. You may install and use the DA module on another computer, but the DA module may not be in use on more than one computer at a time unless you purchase additional licenses. You may make back-up copies of the software for archival purposes. You may permanently transfer your license to use the software to another party who will be bound by this agreement, provided you do not retain any copies of the software.
3. The software is protected by the copyright laws of the U.S. and other countries, and we retain all intellectual property rights in the software. You may not separately publish, sell, market, distribute, lend, lease, rent, or sublicense the software code. However, this license is not to be construed as prohibiting or limiting any fair use sanctioned by copyright law, such as permitted library and classroom usage.
4. We warrant that the software will provide the features and functions generally described in the product specification on our web site when you purchased it and in the product documentation. Media on which the Software is furnished, if any, will be free from defects in materials and workmanship.
5. We have taken all reasonable steps to keep the software free of viruses, spy ware, "back door" entrances, or any other harmful code. We will not track or collect any information about you, your data, or your use of the software except as you specifically authorize. The software will not download or install patches, upgrades, or any third party software without getting your permission. We will not intentionally deprive you of your access to your data.
6. We do not warrant that the software or your ability to use it will be uninterrupted or error-free. To the extent permitted by applicable law, we disclaim any implied warranty of merchantability or fitness for a particular purpose.
7. Your exclusive remedy under the above limited warranty shall be, at our option, either a full refund of the purchase price or correction of the defective software or media. To the fullest extent permitted by applicable law, we disclaim all liability for indirect or consequential damages that arise under this license agreement. Nothing in this agreement limits our liability to you in the event of death or personal injury resulting from gross negligence, fraud, or knowing misrepresentation on our part.
8. If any part of this agreement is found to be invalid or unenforceable, the remaining terms will stay in effect. This agreement does not prejudice the statutory rights of any party dealing as a consumer.
9. This agreement is governed by the laws, including Article 2 of the Uniform Commercial Code, of Massachusetts.
10. This agreement does not supersede any express warranties we made to you. Any modification to this agreement must be agreed to in writing by both parties.

www.diagnosisilk.com

US: Diagnosis LLC; 55 Technology Drive, Suite 100, Lowell, MA 01851; 978-458-1600; sales@diagnosisilk.com
EU: Diagnosis Vision Ltd; Office 117, DOC Building, Balheary Road, Swords, Dublin, K67 6SA0, Ireland; +44 (0) 1223 520699; mail@diagnosisvision.com
UK: Diagnosis UK Ltd; 5 Trust Court, Chivers Way, Vision Park Histon, Cambridge, CB24 9PW, UK; +44 (0) 1223 520699; mail@diagnosisuk.co.uk

The information listed herein is accurate as of the date of printing, however may change at any time without notice. The contents may differ from the current status of approval of the product in your country. Please contact your local Diagnosis representative for more information. © 2021 by Diagnosis, LLC. All rights reserved.
Doc: 14723 Rev: D ECH 1805 Date: 15 Feb 2021

Test Report



Clear Eye Clinic, Lowell, MA

Patient: JEFF FARMER

DOB: 1/11/1968

Age: 52 Years


Gender: Male

Test date: 5/7/2020

Right Eye (OD)

Results Table

# Eye	Color	Date	Results	Rejects	EB%
1 RE	Red	5/7/2020	31	0	0
2 RE	Green	5/7/2020	31	1	0
# Eye	Threshold	Log cd s/m ²	Time	Type	Mode
7 RE [CT]	Cone Threshold	-3.44	11:39	Absolute	Manual
8 RE [RT]	Rod Threshold	-6.91	27:15	Absolute	Manual

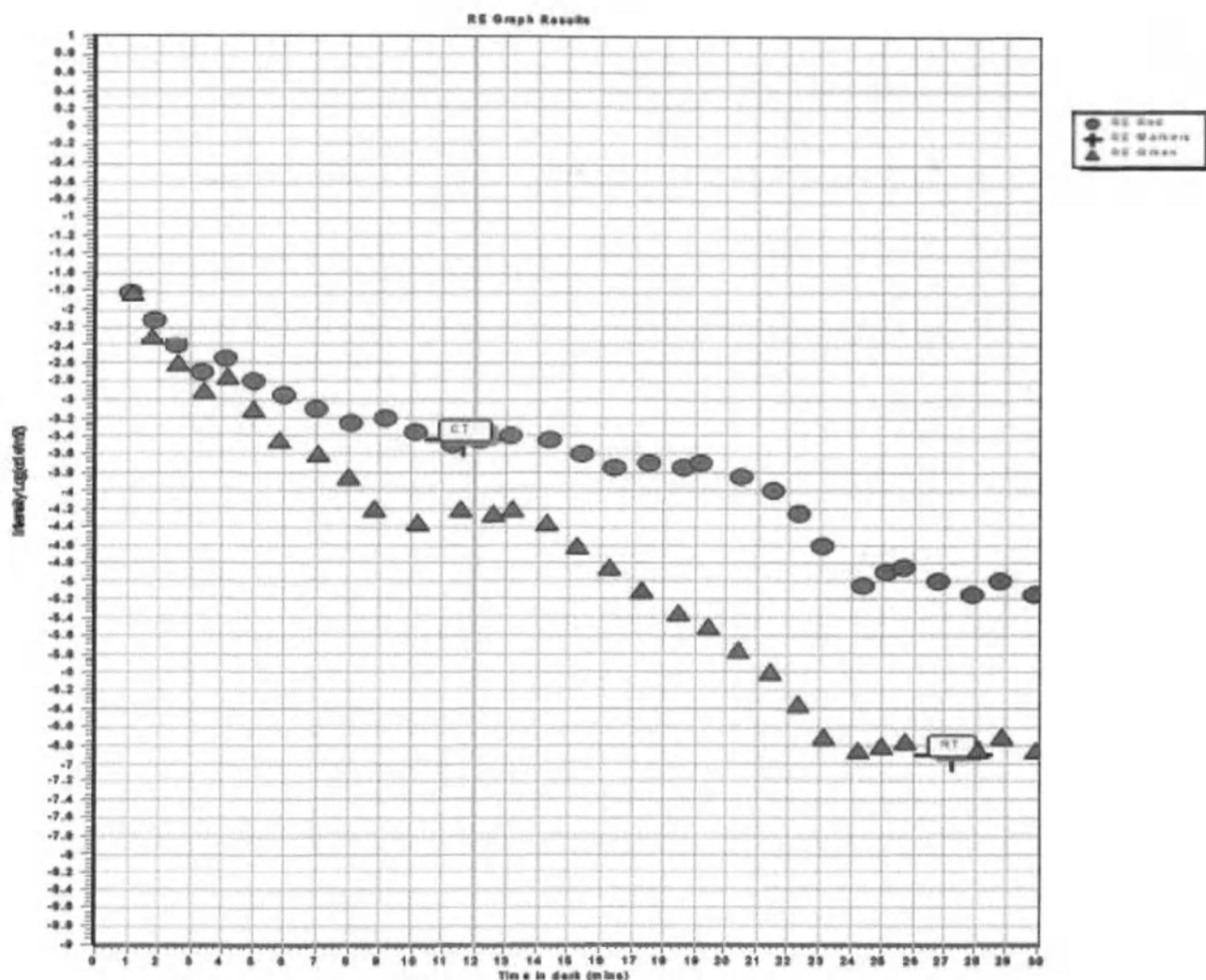

Test comment: jeff full test; full field RE; full bleach

Test type: Cone and Rod Dark AdaptometryC RevC 15515-B ECN 1395 28 September Field of view:

Pupil diameter (CIC08) mm: 5.0/5.0


Electrodes:


Software GW: 6.83.19825.24, Envoy 9.1N: 0-000


Examiner:

Printed: 5/7/2020 1:14:10 PM

Page: 1 of 1

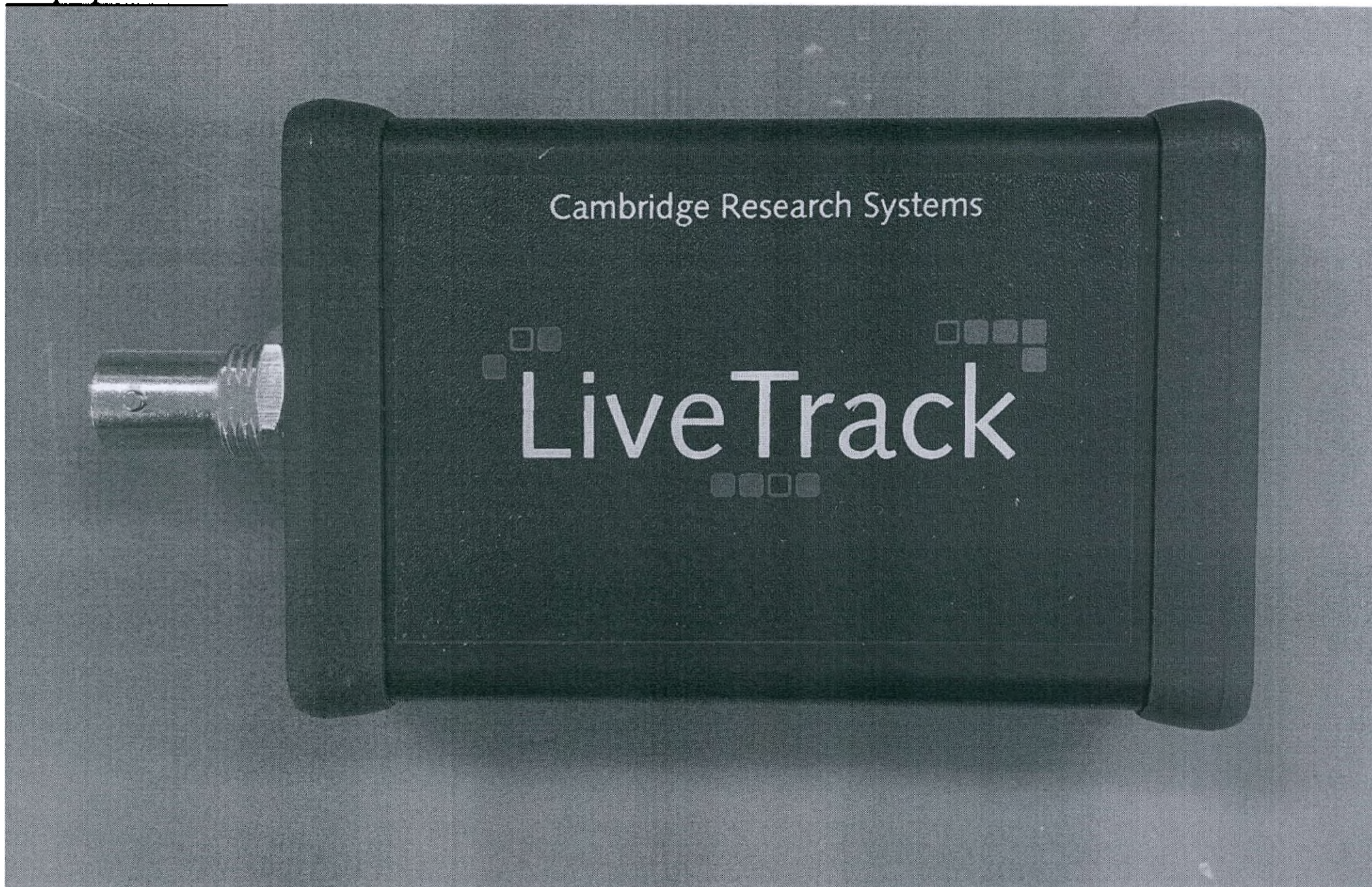

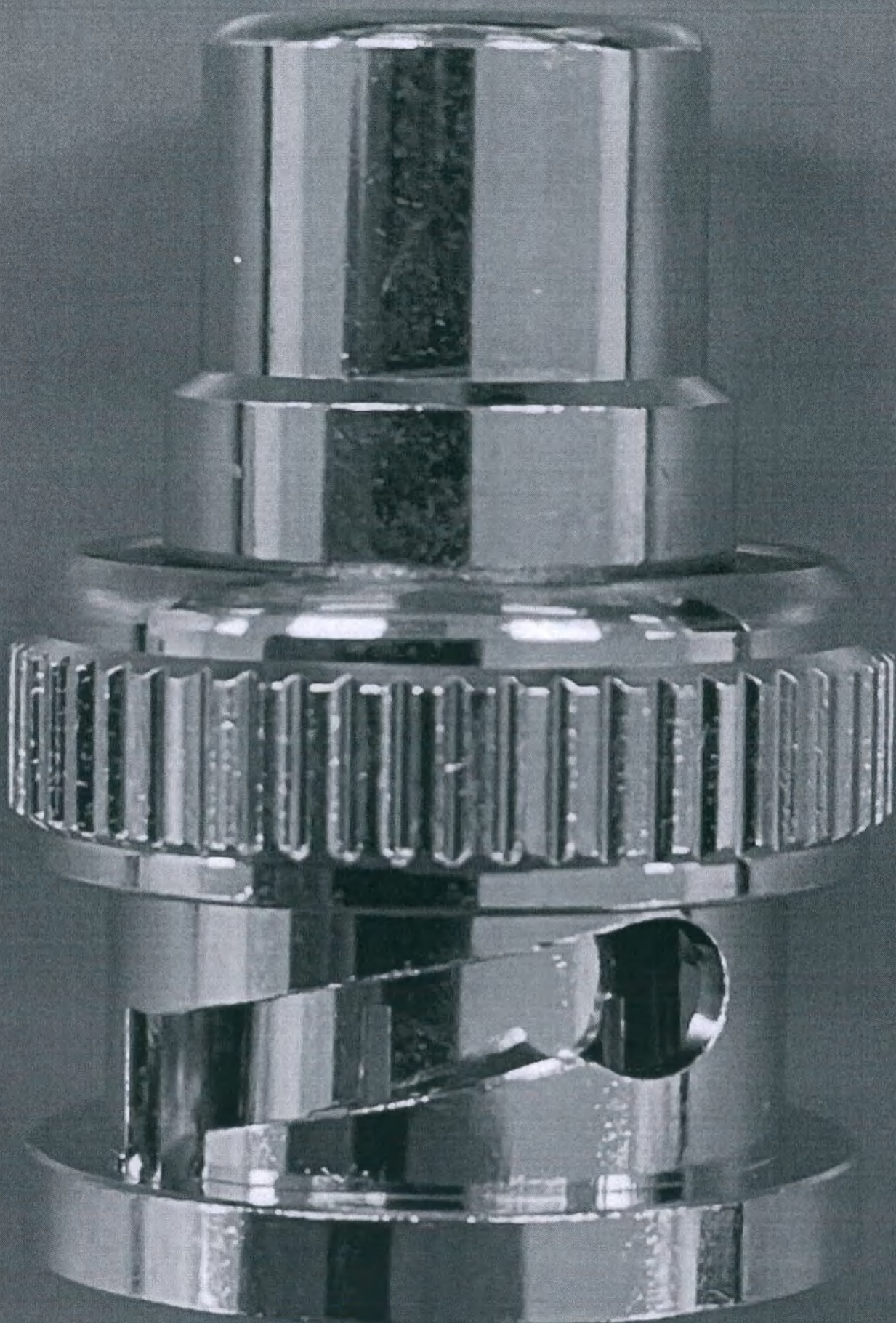

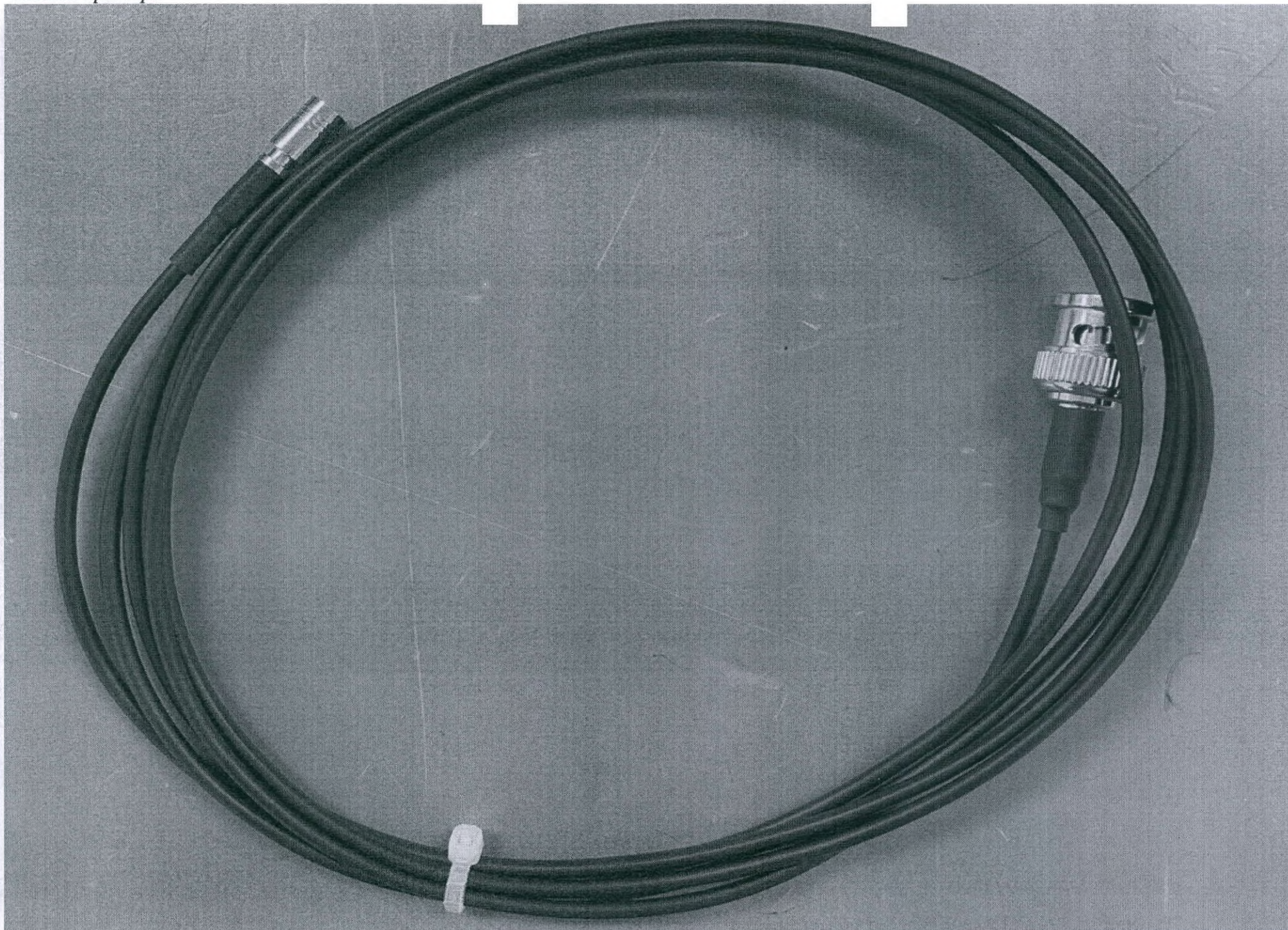



7mm



6mm



5mm



4mm



3mm



2mm




EYETRACKER
TOOLBOX

Scale calibration stick



CAMBRIDGE RESEARCH SYSTEMS

- USB-кабель

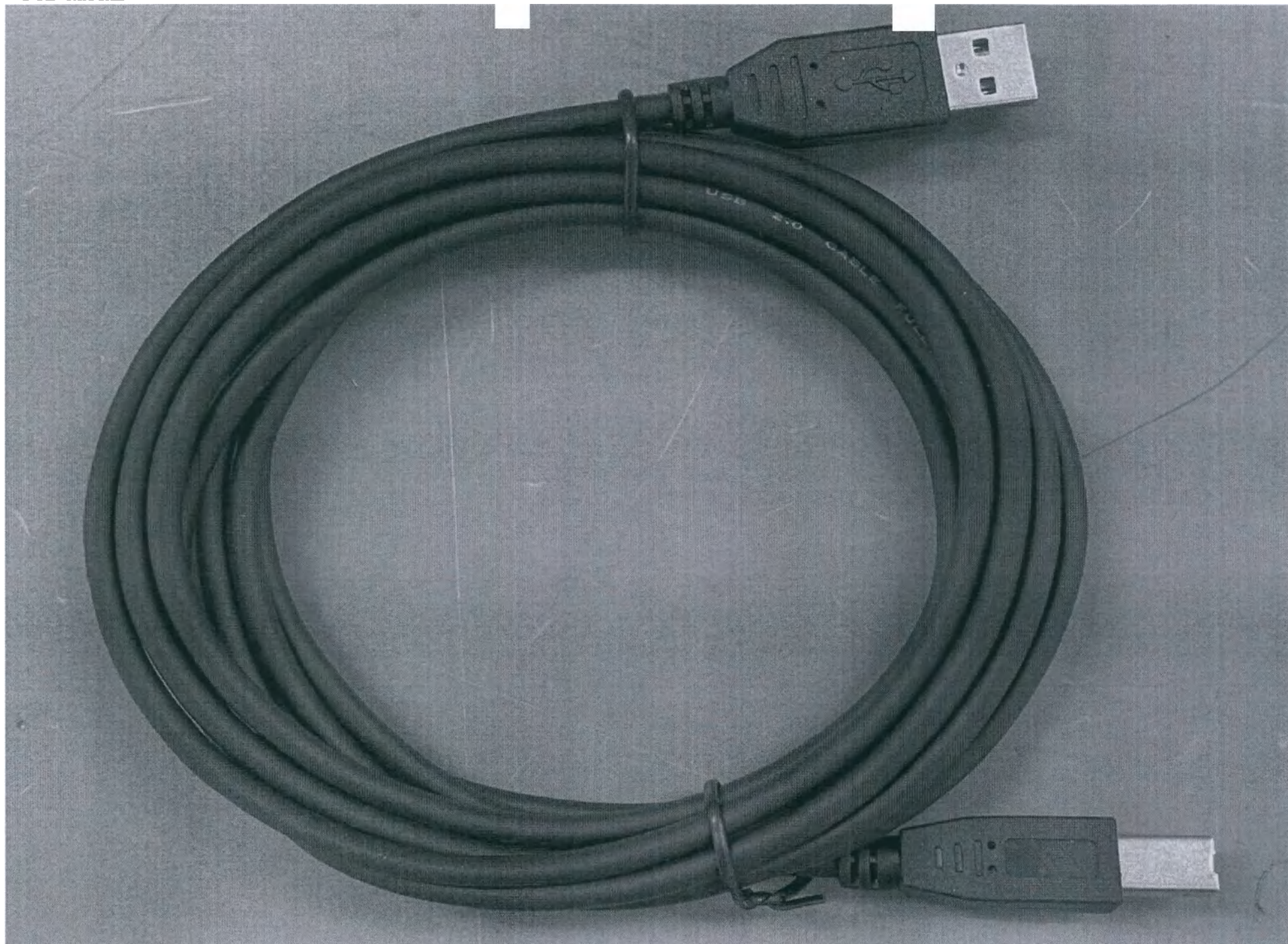




Patient: DOB:12/30/1999 Male

Examiner: SYSDBA Date Tested: 8/30/2018 (Updated: 1/21/2019) Time Tested: 10:48:40 AM


Test Comment: naohiro matozera Unique ID: 37E3A52-9EB9-4D85-A879-2EED-CA07-7D72

Pupillometer.Dia.WRB.BOTH Eyes.Rev H


1 - 1000ms 0.001 (P)cd/m² Blue 0 (P)cd/m² White-6500K

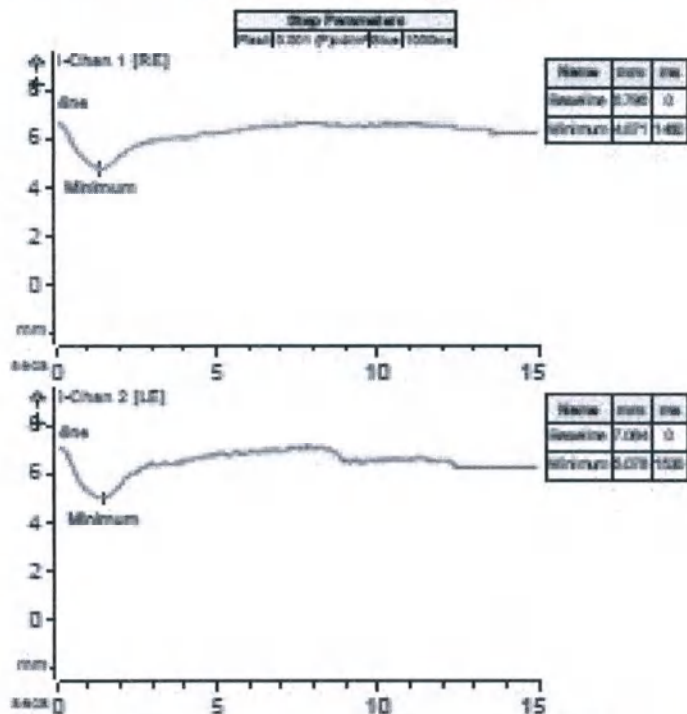


2 - 1000ms 10 (P)cd/m² Red 0.1 (P)cd/m² Blue

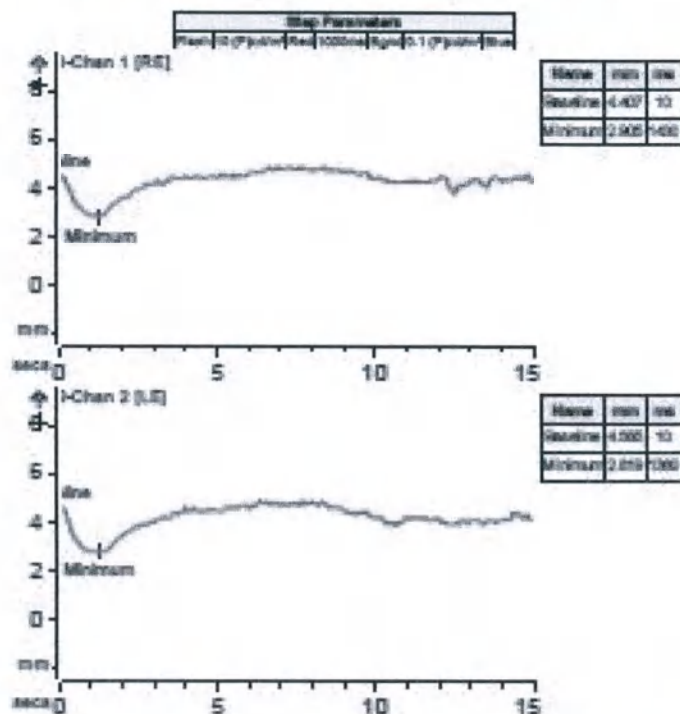


3 - 4ms 3 (P)cd/m² White-6500K 1 (P)cd/m² White-6500K

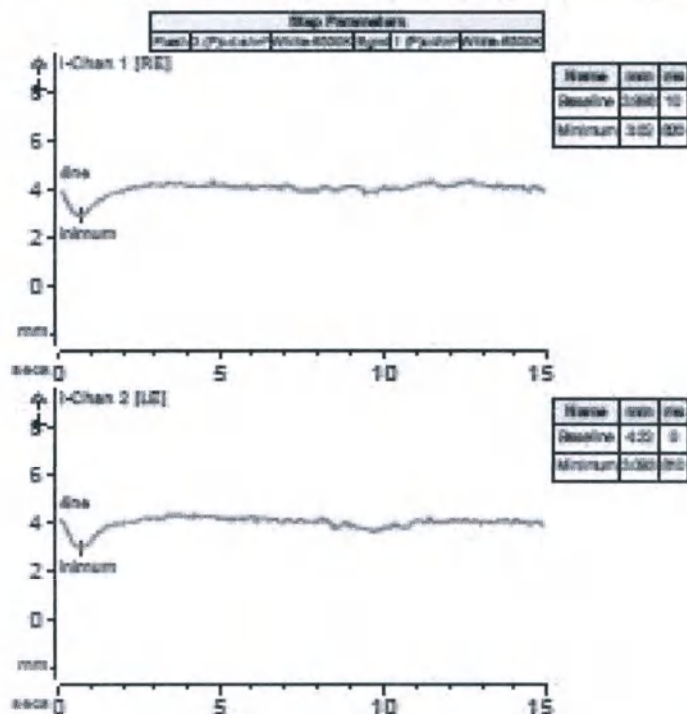


4 - 1000ms 150 (P)cd/m² Blue 0 (P)cd/m² White-6500K

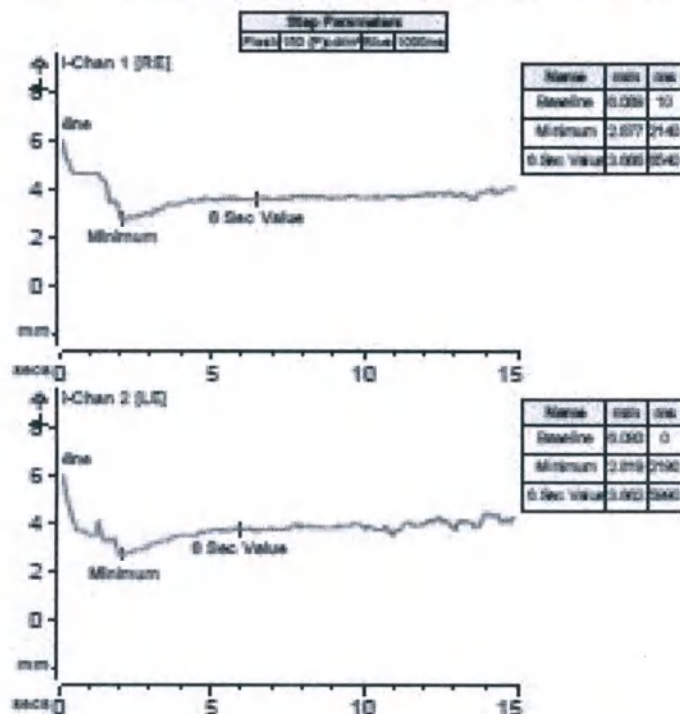

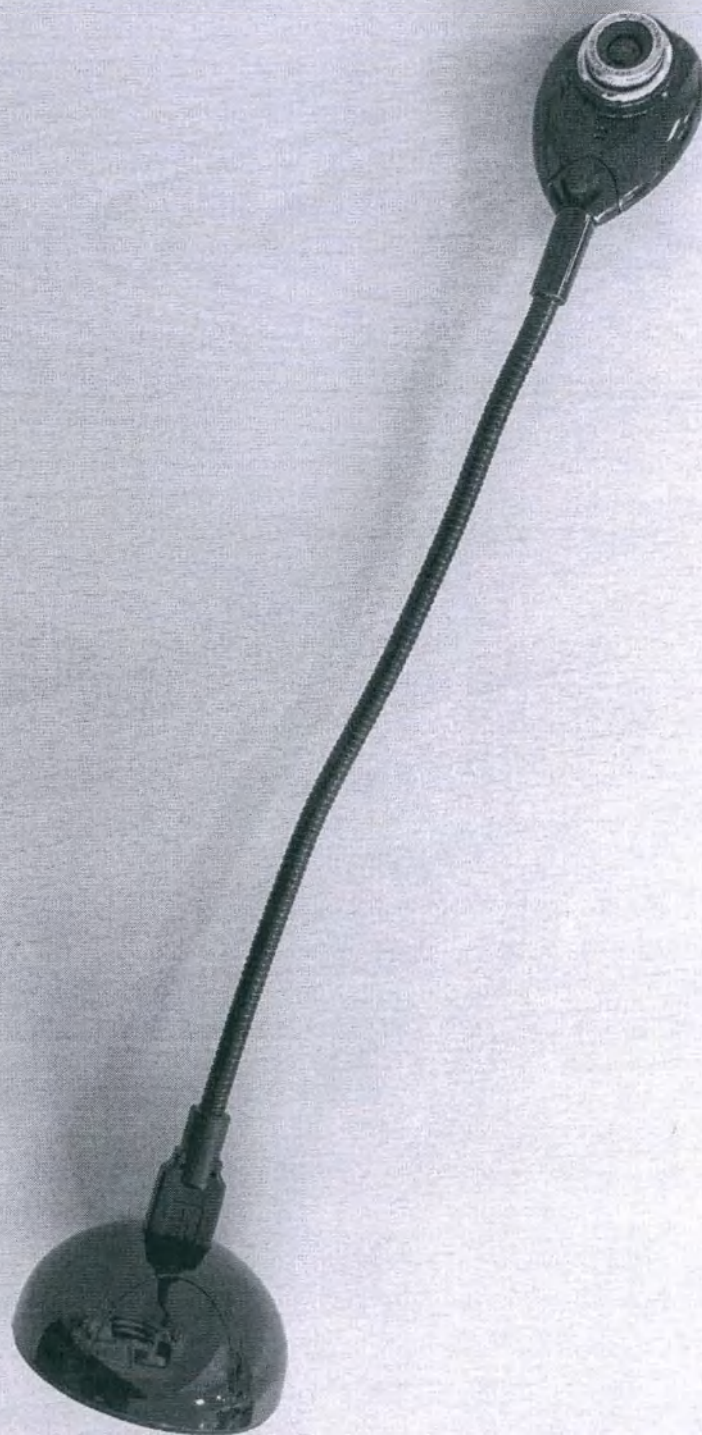





- Слайды с часто используемыми полями зрения

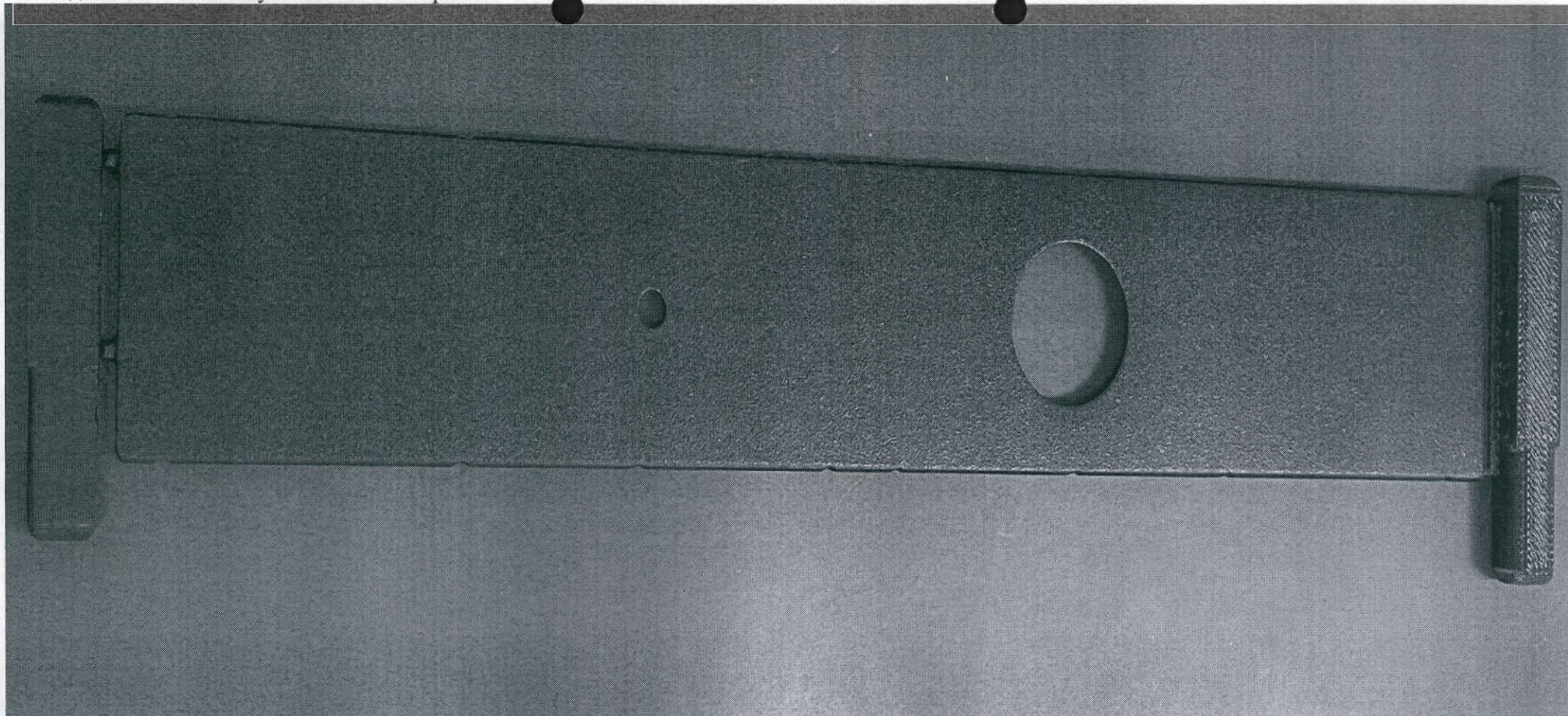

- Кабель для подключения iMask к стимулятору ColorDome

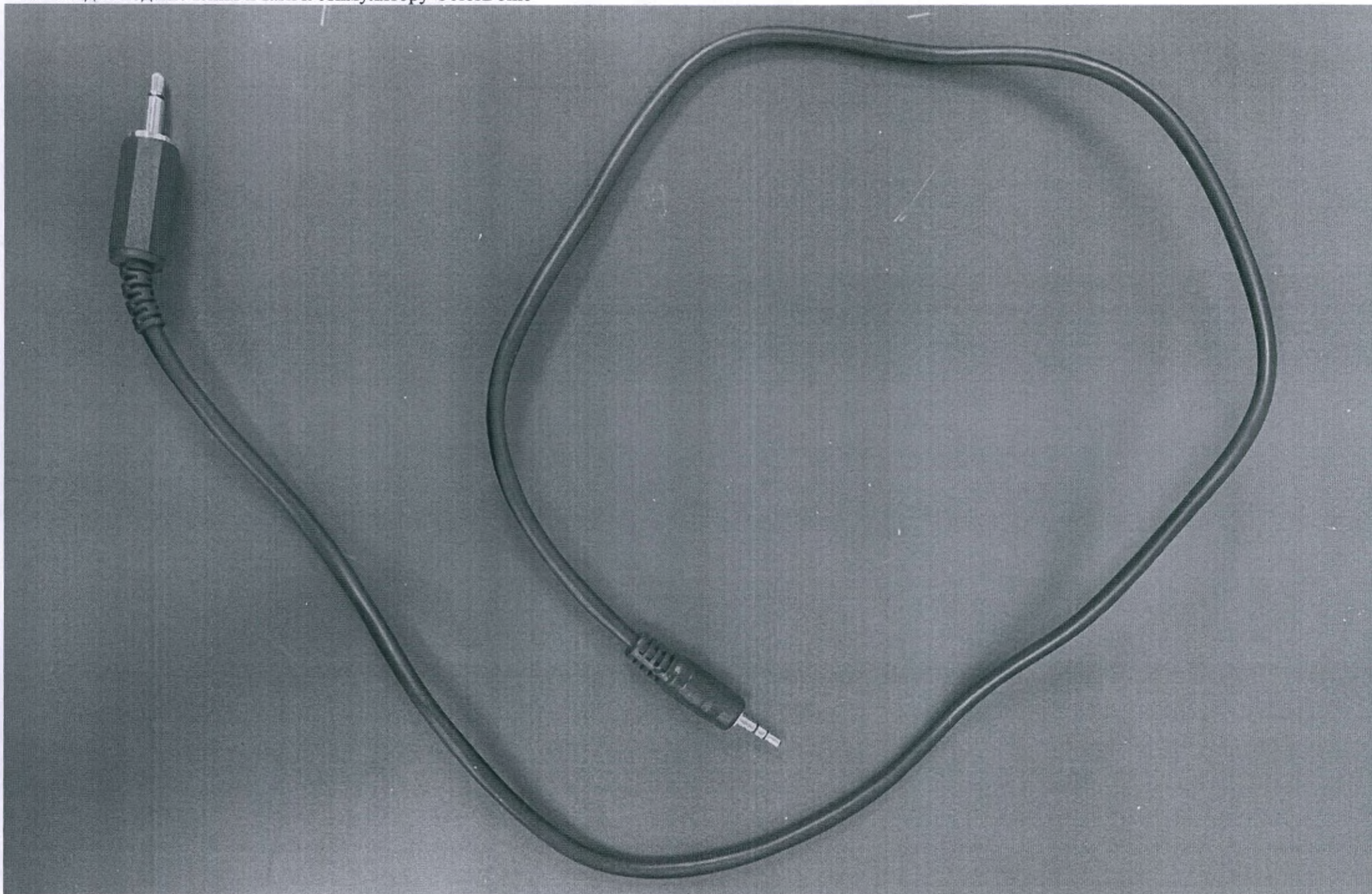

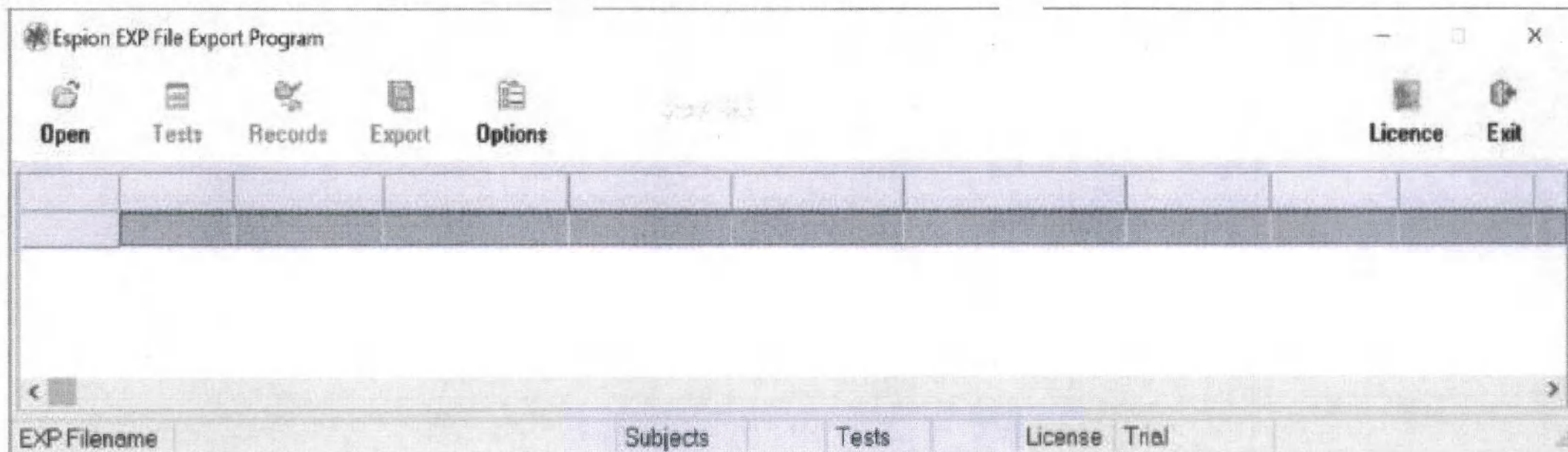

Sweep VEP Module: License Key



Sweep VEP Module License

License Key: XXXX-XXXXXX

License holder:

Name: _____

Address: _____

Date: _____

Fair Use License & Limited Warranty

1. We grant you one license to install and use this software on a single computer at a time. If you do not agree to the following terms of this license, please uninstall and remove all copies and return the product to the place that you purchased it from within 30 days of your purchase for a full refund.
2. You may install and use the module on another computer, but the module may not be in use on more than one computer at a time unless you purchase additional licenses. You may make back-up copies of the software for archival purposes. You may permanently transfer your license to use the software to another party who will be bound by this agreement, provided you do not retain any copies of the software.
3. The software is protected by the copyright laws of the U.S. and other countries, and we retain all intellectual property rights in the software. You may not separately publish, sell, market, distribute, lend, lease, rent, or sublicense the software code. However, this license is not to be construed as prohibiting or limiting any fair use sanctioned by copyright law, such as permitted library and classroom usage.
4. We warrant that the software will provide the features and functions generally described in the product specification on our web site when you purchased it and in the product documentation. Media on which the Software is furnished, if any, will be free from defects in materials and workmanship.
5. We have taken all reasonable steps to keep the software free of viruses, spy ware, "back door" entrances, or any other harmful code. We will not track or collect any information about you, your data, or your use of the software except as you specifically authorize. The software will not download or install patches, upgrades, or any third party software without getting your permission. We will not intentionally deprive you of your access to your data.
6. We do not warrant that the software or your ability to use it will be uninterrupted or error-free. To the extent permitted by applicable law, we disclaim any implied warranty of merchantability or fitness for a particular purpose.
7. Your exclusive remedy under the above limited warranty shall be, at our option, either a full refund of the purchase price or correction of the defective software or media. To the fullest extent permitted by applicable law, we disclaim all liability for indirect or consequential damages that arise under this license agreement. Nothing in this agreement limits our liability to you in the event of death or personal injury resulting from gross negligence, fraud, or knowing misrepresentation on our part.
8. If any part of this agreement is found to be invalid or unenforceable, the remaining terms will stay in effect. This agreement does not prejudice the statutory rights of any party dealing as a consumer.
9. This agreement is governed by the laws, including Article 2 of the Uniform Commercial Code, of Massachusetts.
10. This agreement does not supersede any express warranties we made to you. Any modification to this agreement must be agreed to in writing by both parties.

www.diagnosysllc.com

US: Diagnosys LLC; 55 Technology Drive, Suite 100, Lowell, MA 01851; 978-458-1600; sales@diagnosysllc.com

EU: Diagnosys Vision Ltd; Office 117, DOC Building, Balheary Road, Swords, Dublin, K67 E5A0, Ireland; +44 (0) 1223 520699; mail@diagnosysvision.com

UK: Diagnosys UK Ltd; 5 Trust Court, Chivers Way, Vision Park Histon, Cambridge, CB24 9PW, UK; +44 (0) 1223 520699; mail@diagnosysuk.co.uk

The information listed herein is accurate as of the date of printing, however may change at any time without notice. The contents may differ from the current status of approval of the product in your country. Please contact your local Diagnosys representative for more information. © 2021 by Diagnosys, LLC. All rights reserved.

Doc: 14733 Rev: D 05N 1605 Date: 15 Feb 2021

Stepwise Sweep VEP Report



Clear Eye Clinic, Lowell, MA

Patient:

DOB: 1994

Gender:

Hospital #:


Age: 24 Years

Test date: 8/3/2018

Right Eye (OD)

Snellen 20ft	97.5% Conf.	Calculation	Sfo (CPD)	Log(Sfo)
20/22	> 20/45	Auto	15.71	1.20

RE Response Magnitude vs Log(Spatial Frequency)

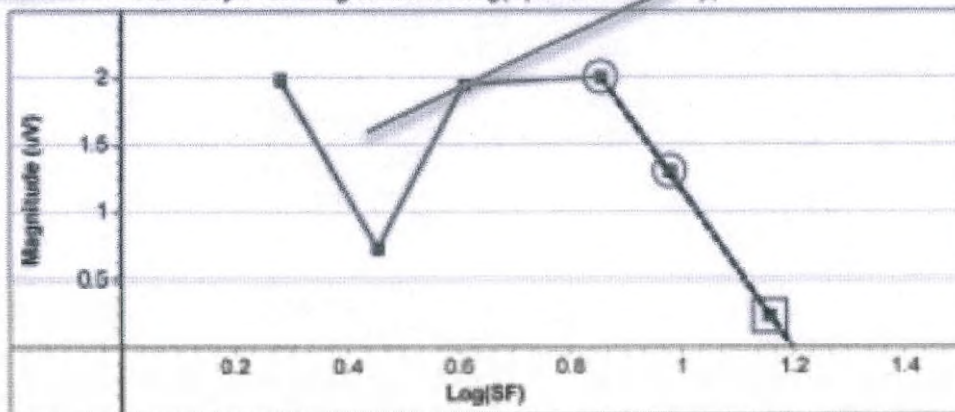


Symbols	
○	Significant
□	Not Significant
◇	Manual

Left Eye (OS)

Snellen 20ft	97.5% Conf.	Calculation	Sfo (CPD)	Log(Sfo)
20/24	> 20/48	Auto	14.66	1.17

LE Response Magnitude vs Log(Spatial Frequency)

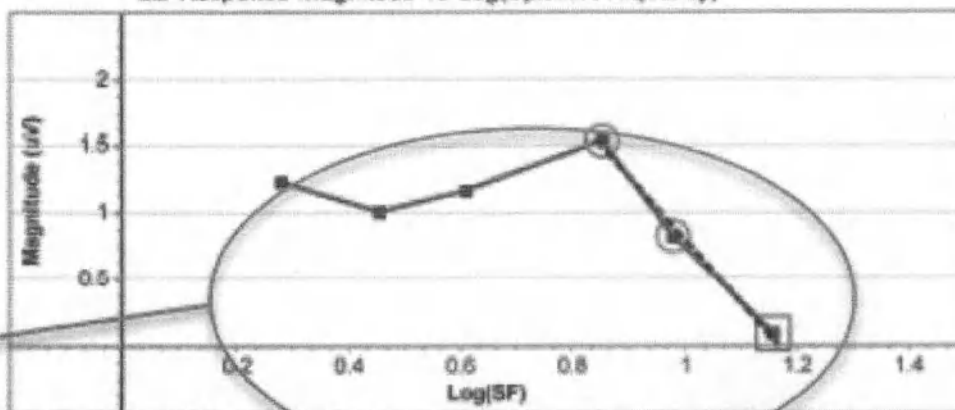

Symbols	
○	Significant
□	Not Significant
◇	Manual

Test comment: bare eye

Test type: Sweep VEP Normal Vision 27m Rev 9.3 (07/17)

Field diameter: 30° (30° mm: 60° 20)

Stimulus: 500 Hz 1000 Hz 1000 Hz 1000 Hz

Field of view: 30° wide x 15° high (diagonal)

Stimulus: Gold Cup Gold Cup


Examiner:

Printed: 5/12/2020 4:07:51 AM

Page: 1 of 1

diagnosys

«На гребне волны»

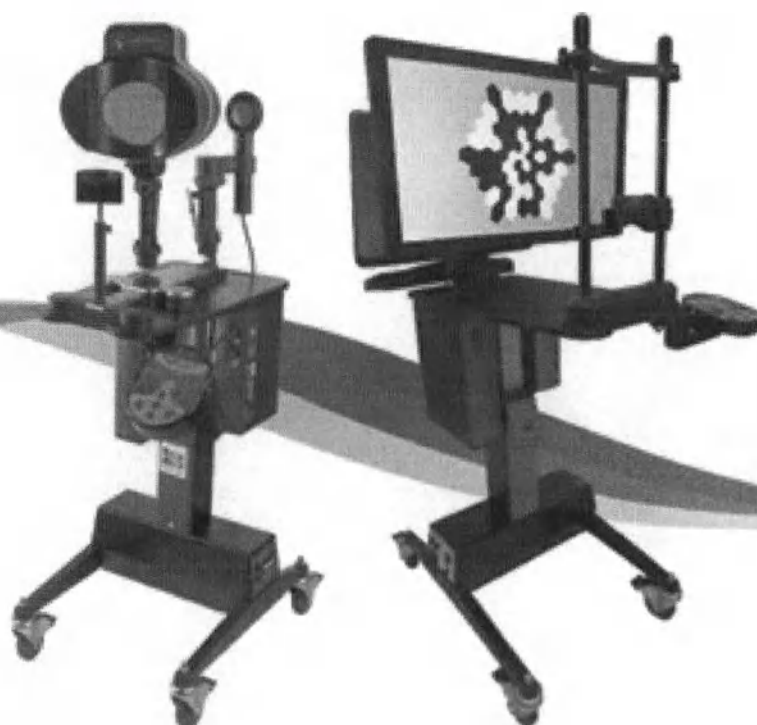

Руководство пользователя Системы Espion Profile



diagnosys LLC, 601114, Nizhny Novgorod region, Nizhny Novgorod, ul. Gagarina, 100/1, 100/1, 100/1

Тел: +7 (831) 418-11-11, 418-11-12, 418-11-13, 418-11-14, 418-11-15, 418-11-16, 418-11-17, 418-11-18, 418-11-19, 418-11-20, 418-11-21, 418-11-22, 418-11-23, 418-11-24, 418-11-25, 418-11-26, 418-11-27, 418-11-28, 418-11-29, 418-11-30, 418-11-31, 418-11-32, 418-11-33, 418-11-34, 418-11-35, 418-11-36, 418-11-37, 418-11-38, 418-11-39, 418-11-40, 418-11-41, 418-11-42, 418-11-43, 418-11-44, 418-11-45, 418-11-46, 418-11-47, 418-11-48, 418-11-49, 418-11-50, 418-11-51, 418-11-52, 418-11-53, 418-11-54, 418-11-55, 418-11-56, 418-11-57, 418-11-58, 418-11-59, 418-11-60, 418-11-61, 418-11-62, 418-11-63, 418-11-64, 418-11-65, 418-11-66, 418-11-67, 418-11-68, 418-11-69, 418-11-70, 418-11-71, 418-11-72, 418-11-73, 418-11-74, 418-11-75, 418-11-76, 418-11-77, 418-11-78, 418-11-79, 418-11-80, 418-11-81, 418-11-82, 418-11-83, 418-11-84, 418-11-85, 418-11-86, 418-11-87, 418-11-88, 418-11-89, 418-11-90, 418-11-91, 418-11-92, 418-11-93, 418-11-94, 418-11-95, 418-11-96, 418-11-97, 418-11-98, 418-11-99, 418-11-100

www.espiondiagnosys.com



colorburst™



pattern



colordome™



envoy





colorflash

espion^e


ДОКУМЕНТ № 15685 | РЕД.С | от 12/8/2020 | БОНА 1565

Электрод DTL Plus – не более 50 шт.

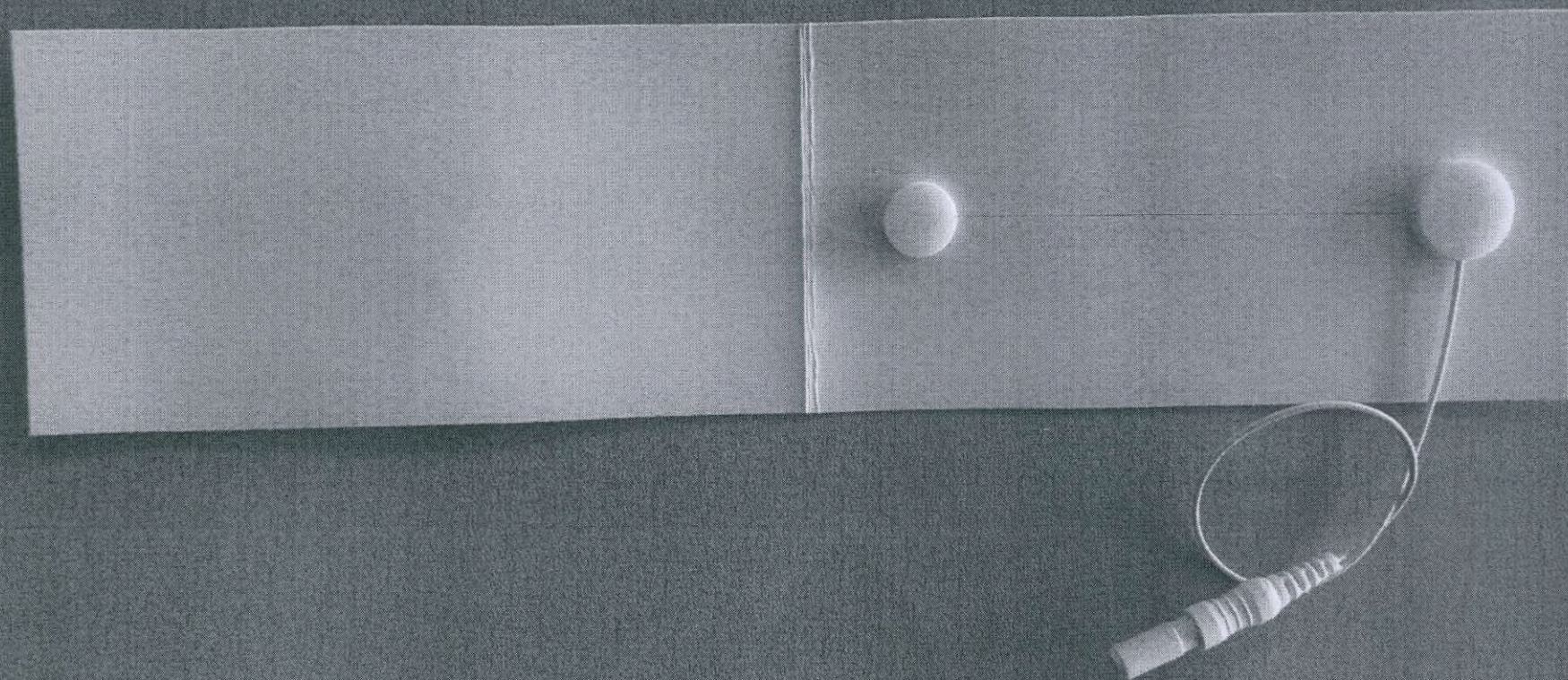

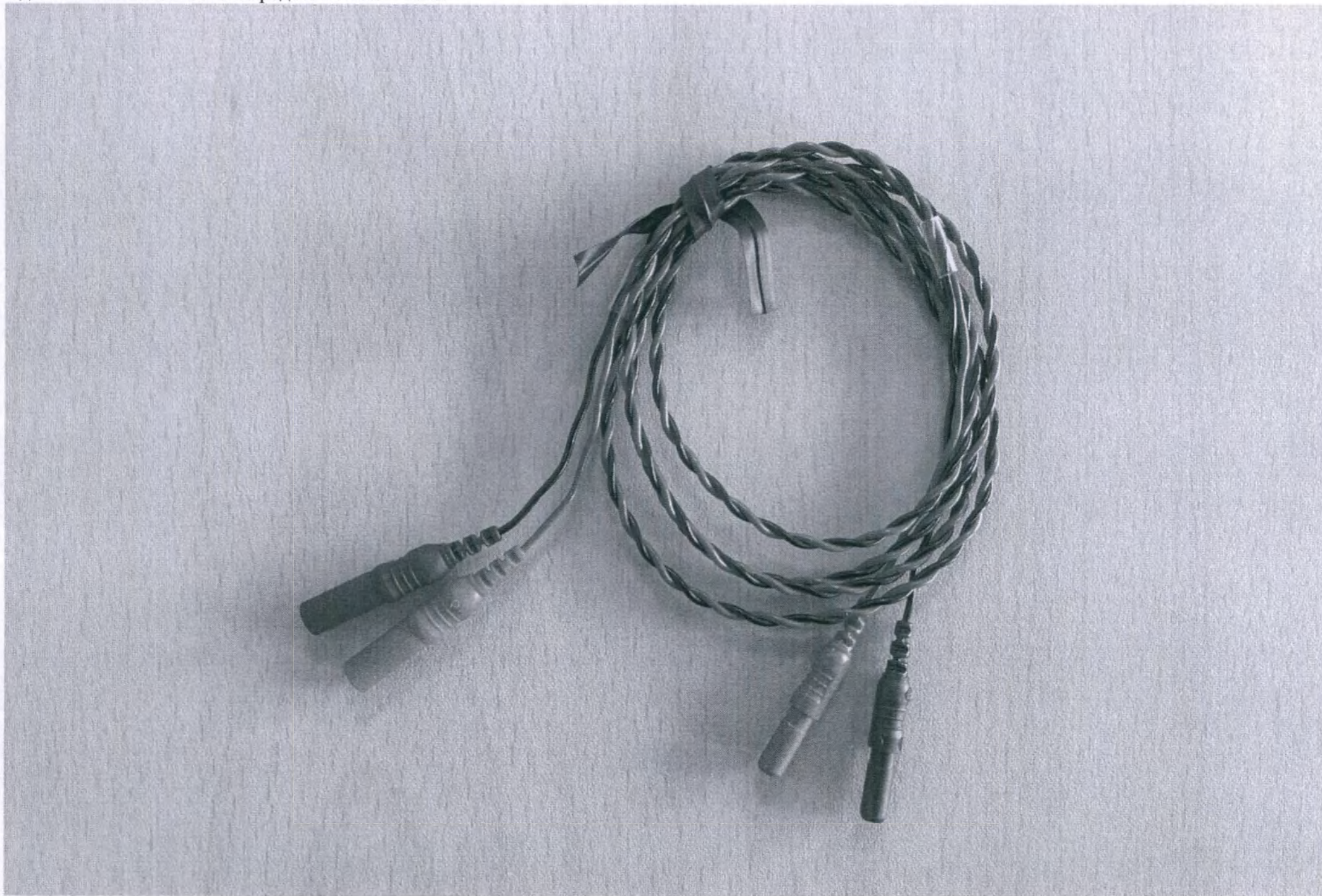

Красный фильтр (22" x 17")

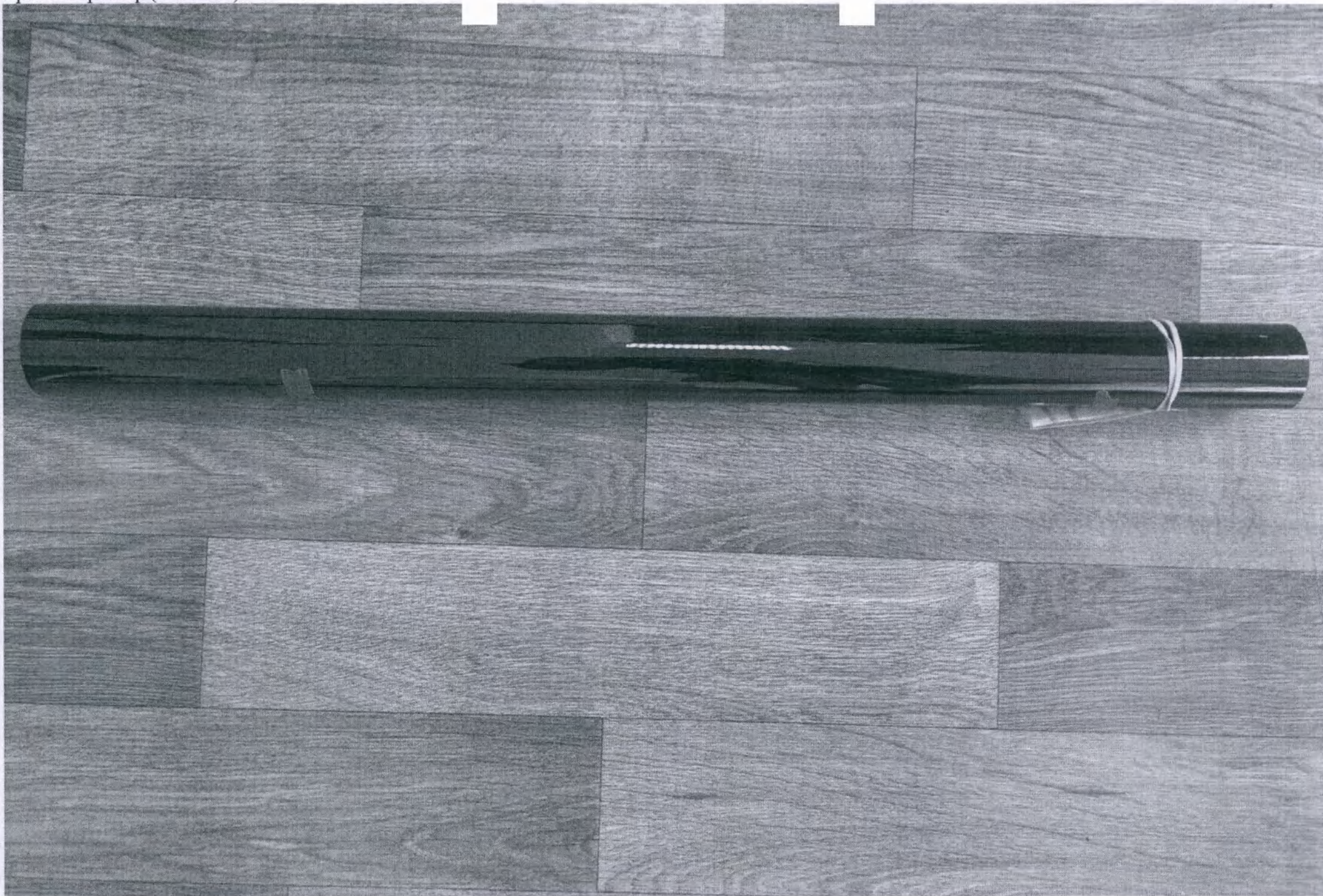









Гаечный ключ – 2 шт.



Набор складных ключей со сферическим концом – 5 шт./набор

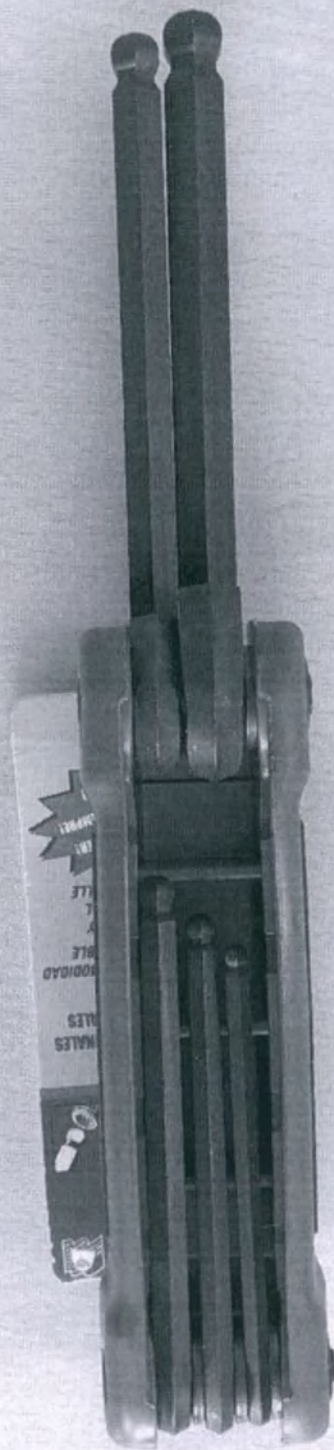




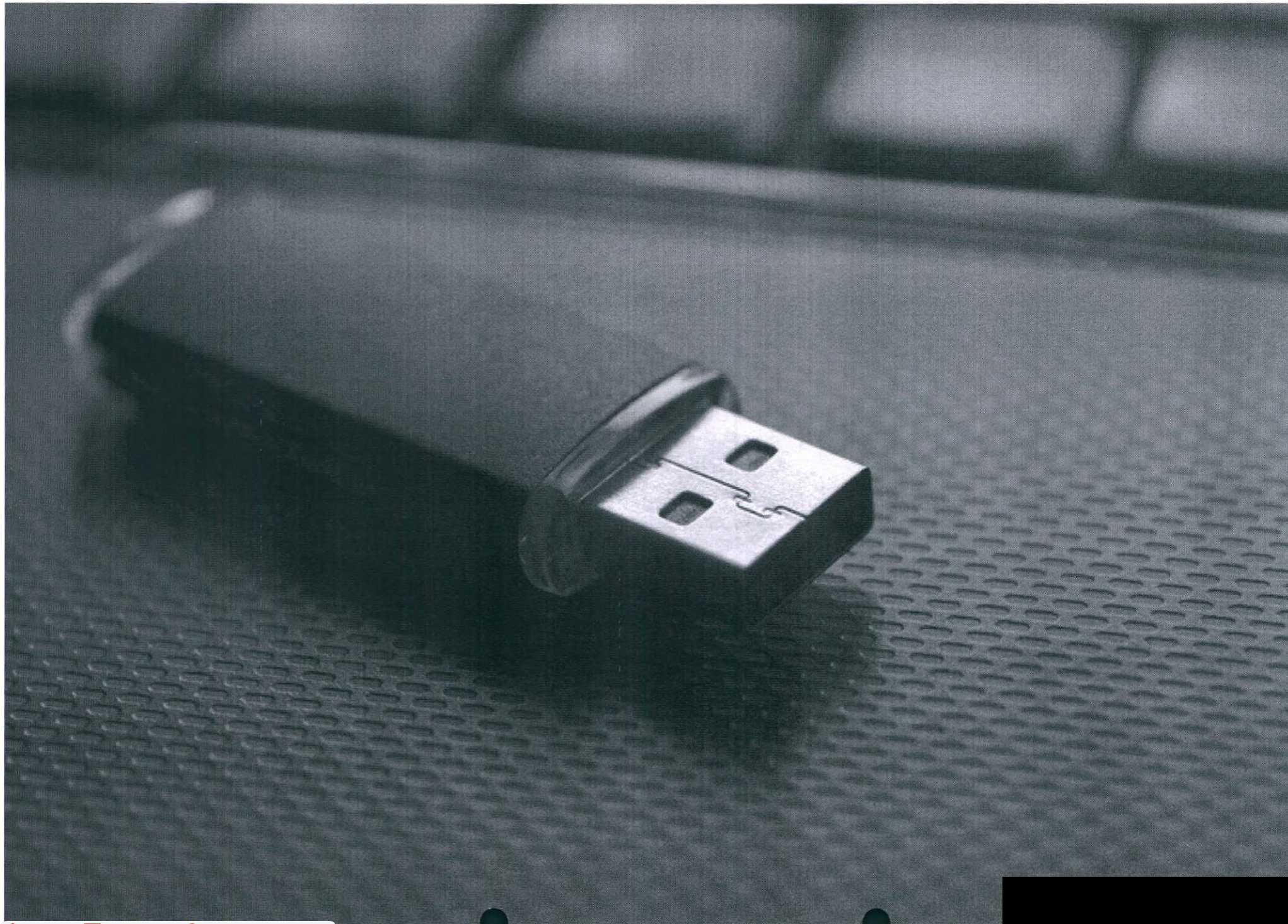




Общий вид системы Espion Profile в конфигурации Blue



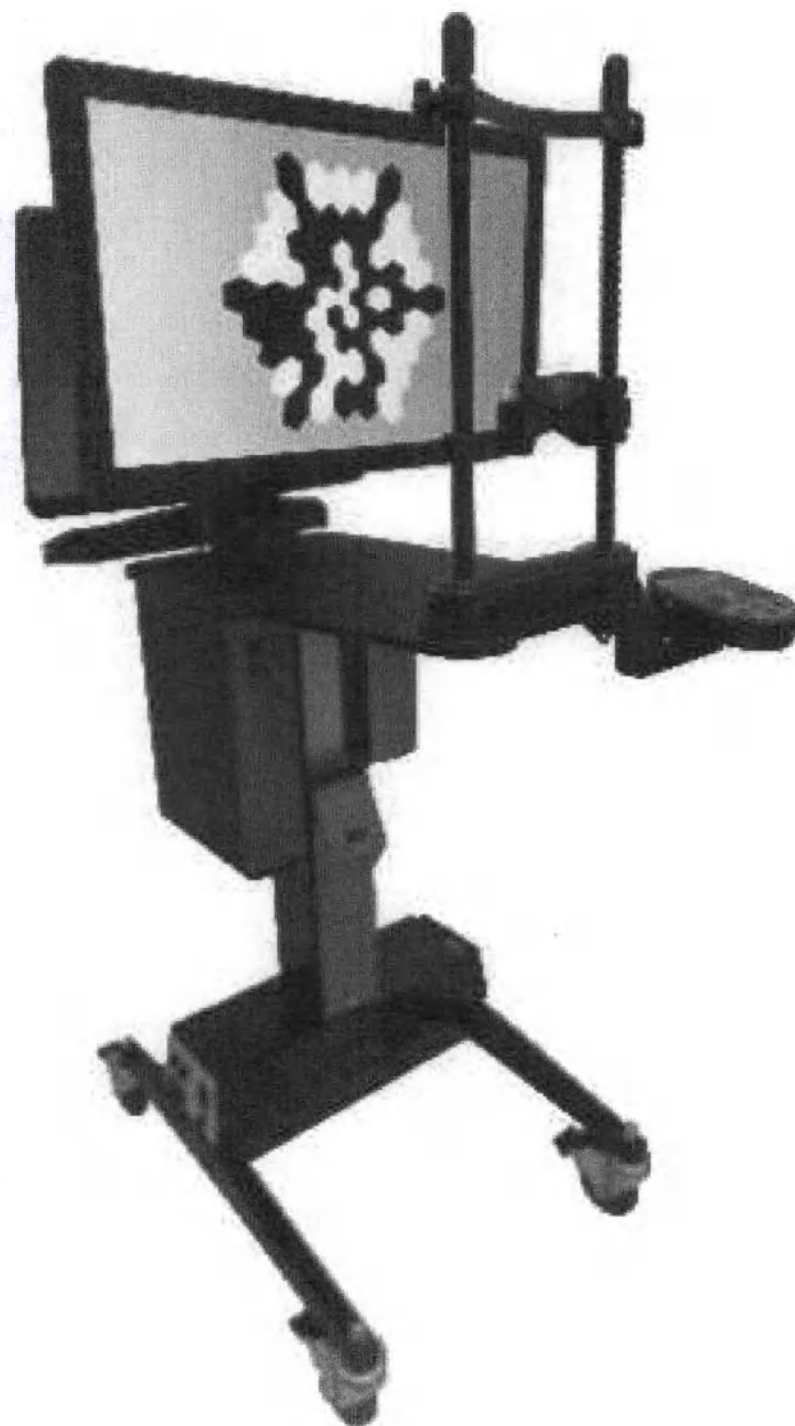

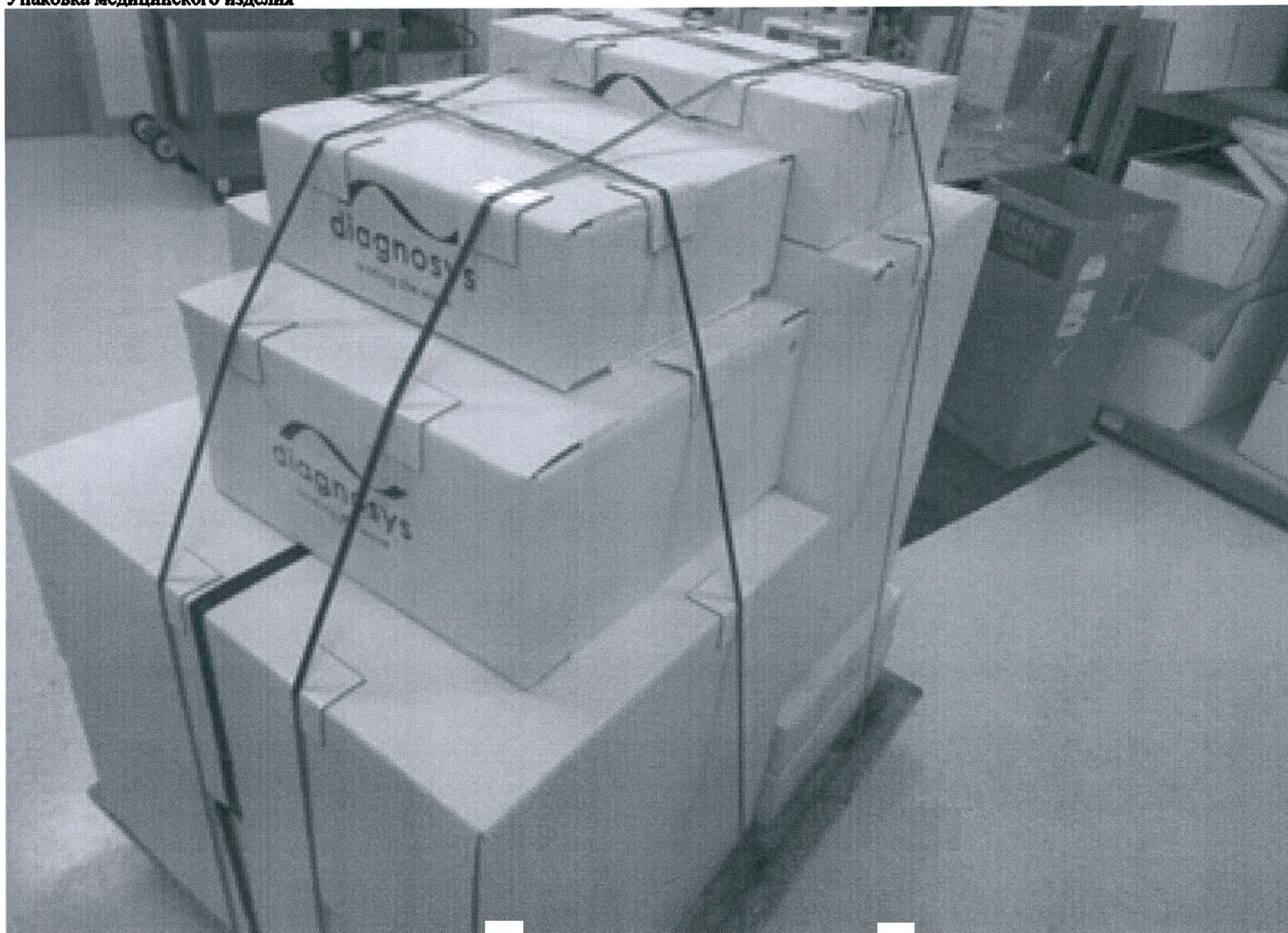



Handwritten signature in blue ink.

Handwritten text in blue ink, likely a date or reference number.