**Medtronic***Alleviating Pain · Restoring Health · Extending Life***CERTIFICATE OF AUTHENTICITY
AND SIGNATURE WITNESS**

State of Massachusetts
County of Essex

On this 18th day of November 2021 before me, the undersigned Notary Public, personally appeared Varsha Kamat, proved to me through satisfactory evidence of identification, which were their driver's license, to be the person whose name is signed on the preceding or attached document in my presence.

(Official signature and seal of Notary)

Signature of Document Custodian (Varsha Kamat)

Medtronic, Inc.
37A Cherry Hill Drive
Danvers, MA 01923

Subscribed and sworn to before me this the
18th day of November 2021

Shannon M. Duffy

My commission expires on April 12, 2024

Medtronic

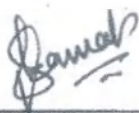
Medtronic, Inc.
710 Medtronic Parkway
Minneapolis, MN 55432
USA
www.medtronic.com

tel 001-763-514-4000
fax 001-763-514-4879

To Whom It May Concern

Medtronic Inc., a corporation organized under the laws of the United States, having its principal place of business at 710 Medtronic Parkway, Minneapolis, MN 55432, USA, as the legal manufacturer of 'DxTerity Diagnostic Angiographic Catheter' hereby confirm that the attached document - Instruction for use M052478T001 Rev. 1D is a true and original document.

Yours sincerely,

Signed on behalf of Medtronic Inc.

Varsha Kamat
Sr. RA Specialist

Nov 18, 2021

Date:

Medtronic

DxTerity™

Diagnostic Catheter
Диагностичен катетър
Diagnostický katéter
Diagnostisk kateter
Diagnostikkatheter
Διαγνωτικός καθετήρας
Catéter de diagnóstico
Diagnostiline kateeter
Diagnostinen katetri
Cathéter de diagnostic
Dijagnostički kateter
Diagnosztikai katéter
Catetere diagnostic
진단용 카테터
Diagnostinis kateteris
Diagnostikas katetrs
Diagnostische katheter
Diagnostisk kateter
Cewnik diagnostyczny
Cateter de diagnóstico
Cateter de diagnóstico
Cateter de diagnostic
Диагностический катетер
Diagnostický katéter
Diagnostični kateter
Dijagnostički kateter
Diagnostisk kateter
Dijagnostik Kateter
Диагностичний катетер
诊断导管

Instructions for Use • Инструкции за употреба •
Návod k použití • Brugsanvisning • Gebrauchsanweisung •
Οδηγίες χρήσης • Instrucciones de uso • Kasutusjuhend •
Käyttöohjeet • Mode d'emploi • Upute za upotrebu •
Használati útmutató • Istruzioni per l'uso • 사용 지침 •
Naudojimo instrukcijos • Lietošanas pamācība •
Gebruiksaanwijzing • Bruksanvisning •
Instrukcja użytkowania • Instruções de utilização •
Instruções de utilização • Instrucțiuni de utilizare •
Инструкция по эксплуатации • Pokyny na používanie •
Navodila za uporabo • Uputstva za upotrebu •
Bruksanvisning • Kullanım Talimatları •
Інструкція з експлуатації • 使用说明

2016MAY/27 at 11:28 a.m. Doc number: M052478T001 Rev. 1D [global]
Printspec v - 3.5 x 10.75 inches

master:20160605, altDocMod:20160504, altDocMod2:20151000, altDocStyle_graph:
20160504, altDocStyle:20160504, master_pageSet:20160504, languageLookups:20160530,
table:20160113, characterJoin:20110211, Join:20110211, Join_rev:20110211

Explanation of symbols on package labeling /
Обяснение на символите от етикетите на
опаковката / Vysvětlení symbolů uvedených na
štítcích balení / Forklaring af symbolerne på
emballagens mærkater / Erläuterung der Symbole auf
dem Verpackungsetikett / Επεξήγηση των συμβόλων
στη σήμανση της συσκευασίας / Explicación de los
símbolos que aparecen en el etiquetado del envase /
Pakendi siltidel olevate sümbolite tähendus /
Pakkausetiketin symbolien selitykset / Explication des
symboles des étiquettes sur l'emballage / Objašnjenje
simbola na ambalažnoj etiketi / A csomagoláson
látható szimbólumok jelentése / Spiegazione dei
simboli presenti sulle etichette della confezione / 포장
레이블의 기호 설명 / Pakuotės etikečių simbolių
reikšmės / Uz iepakojuma etiķetes attēloto simbolu
skaidrojums / Verklaring van de symbolen op de
verpakingslabels / Forklaring av symboler på
pakningsetiketten / Objaśnienie symboli znajdujących
się na etykietach opakowania / Explicação dos
símbolos nas etiquetas da embalagem / Explicação dos
símbolos na documentação da embalagem / Legenda
simbolurilor cuprinse în etichetele de pe ambalaj /
Объяснение символов на этикетках упаковки /
Vysvetlivky k symbolom na označení balenia / Razlaga
simbolov na embalaži / Objašnjenje simbola na
oznakama na pakovanju / Förklaring av symboler på
förpackningsmärkningen / Ambalaj etiketi üzerindeki
sembollerin açıklaması / Пояснения символов на
маркуванні упаковки / 包装标签上的符号说明

Quantity / Количество / Množství / Antal / Menge /
Ποσότητα / Cantidad / Kogus / Määrä / Quantité /
Količina / Mennyiség / Quantità / 수량 / Kiekis /
Daudzums / Antal / Antall / Ilość / Quantidade /
Quantidade / Cantitate / Количество / Množstvo /
Količina / Količina / Antal / Miktar / Кількість / 数量

Open here / Отворете тук / Zde otevřete / Åbnes
her / Hier öffnen / Ανοίξτε εδώ / Abir aquí / Avage
silt / Avaa tästä / Ouvrir ici / Ovdje otvoriti / Itt nyílik /
Aprire qui / 여기를 여십시오 / Atidaryti čia / Åtvært her / Hier
openen / Åbnes her / Otwierać tutaj / Abir aquí / Abir
aquí / Deschideți aici / Открывать здесь / Tu otvoriti /
Odpriete tutaj / Ovdje otvoriti / Otvorať hár / Buradan
açın / Відкривати тут / 由此处打开

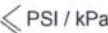
For US audiences only / Само за САЩ / Pouze pro
uživatele z USA / Gælder kun i USA / Gilt nur für
Leser in den USA / Μόνο για πελάτες εντός των
ΗΠΑ / Solo aplicable en EE. UU. / Ainult kasutajatele
Ameerika Ühendriikides / Koskes vain Yhdysvaltoja /
Ne s'applique qu'aux États-Unis / Samo za američko
tržište / Csak egyesült államokbeli felhasználókra /
Eksklusivamente per il mercato statunitense / 미국 내
사용자만 해당 / Tik JAV naudotojams / Tikai pircejiem
ASV / Alleen van toepassing voor de VS / Gælder kun
USA / Dotyczy tylko odbiorców w USA / Aplicable
apenas aos EUA / Aplicável apenas aos EUA / Numai
pentru clienții din SUA / Тільки для США / Len pre
používateľov v USA / Samo za uporabnike v ZDA /
Samo za korisnike u SAD-u / Gæller endast i USA /
Valnizca ABD'deki kullanıcılar için / Лише для
користувачів із США / 仅供美国用户使用

Manufactured in / Произведено в / Vyrobeno v /
Fremstillet i / Hergestellt in / Κατασκευάστηκε σε /
Fabricado em / Tootmiskoht / Valmistuspaikka / Lieu
de fabrication / Zemlja proizvodnje / A gyártás helye /
Fabbriato in / 제조 국가 / Pagaminimo šalis /
Ražošanas vieta / Vervaardigd in / Produziert i /
Miejscę produkcji / Fabricado em / Fabricado em /
Fabricat în / Произведено в / Vyrobené v / Izdelano
v / Proizvedeno u / Tillverkad i / Imalat yeri / Місце
виготовлення / 製造地

Outer diameter / Външен диаметър / Vnější průměr /
Udvendig diameter / Außendurchmesser / Εξωτερική
διάμετρος / Diámetro externo / Välisläbimõõt /
Ulkoläpimitta / Diamètre externe / Vanjski promjer /
Külső átmérő / Diametro esterno / 외부 직경 / Išorinis
skersmuo / Årējais diametrs / Buitendiameter / Ytre
diameter / Średnica zewnętrzna / Diámetro externo /
Diámetro externo / Diametru exterior / Внешний
диаметр / Vonkajši priemer / Zunanji premer / Spoljni
prečnik / Ytterdiameter / Dış çap / Зовнішній
діаметр / 外径

	Usable length / Исползаема дължина / Využitelná délka / Anvendelig længde / Nutzbare Länge / Οφέλιμο μήκος / Longitud útil / Kasutatav pikkus / Käyttöpituus / Longueur utile / Iskoristiva dužina / Hasznos hossz / Lunghezza utile / 사용 가능한 길이 / Naudingas ilgis / Lietderīgais garums / Bruikbare lengte / Nyttbar lengde / Długość użyteczna / Comprimento útil / Comprimento útil / Lungime utilizabilă / Рабочая длина / Použitelná délka / Uporabna dolžina / Upotrebjliva dužina / Arbetstlångd / Kullanılabilir uzunluk / Корисна довжина / 可用长度
	Maximum guidewire diameter / Максимален диаметър на водач / Maximální průměr vodičho drátu / Maksimal guidewirediameter / Maximaler Durchmesser des Führungsdrahts / Μέγιστη διέκτρος οδηγού σύρματος / Diámetro máximo de la guala / Maksimaalne juhteraadi läbimõõt / Ohjainvajerin enimmäisläpimitta / Diamètre maximum du fil-guide / Maksimální průměr žice vodiče / A vezetődírt legnagyobb átmérője / Diametro massimo del filo guida / 가이드 와이어 최대 직경 / Didžiausias kreipiamosios vielos skersmuo / Maksimālais vadītājstīgas diametrs / Maximale diameter voerdraad / Sterste diameter på ledevaier / Maksymalna średnica przewodnika / Diámetro máximo do fio-guia / Diámetro máximo do fio-guia / Diametru maxim fir de ghidare / Максимальный диаметр проводника / Maximálny priemer vodiaceho drôtu / Največji premer vodilne žice / Maksimalni prečnik vodiča-žice / Ledarens maximala diameter / Maksimum kalavuz tel çapı / Максимальный диаметр проводника / 最大导线直径
	Number of sideholes / Брой странични отвори / Počet postranních otvorů / Antal sidehuller / Anzahl der seitlichen Öffnungen / Αριθμός πλευρικών οπών / Número de orificios laterales / Külavade arv / Sivureikien määrä / Nombre d'orifices latéraux / Broj bočnih otvora / Oldalnyílások száma / Numero di fori laterali / 측공 개수 / Soninių angų skaičius / Šonu caurumu skaits / Aantal zijopeningen / Antall sidehull / Liczba otworów bocznych / Número de orificios laterais / Număr de orificii laterale / Количество боковых отверстий / Počet bočných otvorov / Število stranskih odprtín / Broj bočnih otvora / Antal sidohål / Yan deliklerin sayısı / Кількість бокових отворів / 側孔數
	Maximum PSI/kPa / Максимум PSI/kPa / Maximální PSI/kPa / Maksimal PSI/kPa / Maximale psi/kPa / Μέγιστη τιμή PSI/kPa / Valor máximo de presión en PSI/kPa / Maksimaalne PSI/kPa / PSI-kPa-enimmäisarvo / PSI/kPa maximum / Maksimalan PSI/kPa / Maximális PSI/kPa / Pressione massima (PSI/kPa) / 최대 PSI/kPa / Didžiausias slėgis PSI / kPa / Maksimālā PSI/kPa vērtība / Maximale psi/kPa / Maks. psi/kPa / Maksymalne ciśnienie w PSI/kPa / PSI/kPa máximo / PSI/kPa máxima / Valoarea PSI/kPa maximă / Максимально фунтов на кв. дюйм/кПа / Maximálny tlak psi/kPa / Največja vrednost PSI/kPa / Maksimalni pritisk u PSI/kPa / Högsta PSI/kPa / Maksimum PSI/kPa / Максимальный тиск (psi/кПа) / 最大 PSI/kPa
	Do not bend / Не изъгайте / Neohýbajte / Må ikke bøyes / Nicht knicken / Μην το κάμπτετε / No doblar / Mitte painutada / Älä taivuta / Ne pas plier / Ne savijati / Tilos meghajlítani / Non piegare / 구부리지 마십시오 / Nelenkti / Nelociti / Niet buigen / Skal ikke bøyes / Nie wyginać / Não dobrar / Não dobrar / A nu se îndoi / Не сгибать / Neohýbajte / Ne upogibati / Ne savijati / Får inte böjas / Bükmeysin / Не згинати / 勿折

DxTerity™

Diagnostic Catheter

1. Device description

DxTerity™ diagnostic catheters are available in various configurations and French sizes.

2. Indications

The diagnostic catheter is indicated for cardiac and vascular procedures. It is designed to deliver radiopaque media, guidewires, and therapeutic agents to selected sites in the vascular system. The different configurations of the diagnostic catheter are designed to be used in arteries from access sites such as the radial, brachial, and femoral arteries.

3. Contraindications

None specific to this device.

4. Warnings and precautions

Caution: Read all instructions carefully. Failure to follow the instructions, warnings, and precautions may lead to serious consequences or injury to the patient.

4.1. General

- For single-patient use only. Do not reuse, reprocess, or resterilize this product. Reuse, reprocessing, or resterilization may compromise the structural integrity of the device and/or create a new risk of contamination of the device, which could result in patient injury, illness, or death. Cleaning, disinfection, and resterilization may compromise the essential material and design characteristics of the device and lead to device failure.
- Do not use the diagnostic catheter if it has been damaged or any other anomaly is observed.
- Do not use after the Use-by date.
- Use aseptic technique during the entire procedure.
- **Warning:** Administer appropriate anticoagulation therapy to reduce potential thrombosis. If the patient is not appropriately anticoagulated, thrombus formation may occur. Testing without anticoagulation shows variable amounts of thrombus formation on the device surface.
- After use, this product may be a biohazard. Handle and dispose of all such devices in accordance with accepted medical practice and applicable local, state, and federal laws and regulations.
- Only physicians trained in intravascular angiographic procedures should use or supervise the use of the diagnostic catheter.
- Carefully inspect the diagnostic catheter before use to verify that the size, shape, and condition of the catheter are suitable for the intended procedures.
- Store diagnostic catheters in a cool, dry, and dark area.
- Do not expose diagnostic catheters to UV light.
- The diagnostic catheter is supplied sterile and nonpyrogenic. Do not use the catheter if the catheter package has been opened or damaged before use.
- Do not heat or bend the diagnostic catheter before use.
- Select a diagnostic catheter with optimal tip shape and optimal size, taking into account the access site, the target site, and the patient's anatomy.

4.2. Interactions

- When using a drug or device with the diagnostic catheter, be aware of the properties and characteristics of the drug or device and exercise due caution to avoid damaging the diagnostic catheter.
- Avoid using alcohol, antiseptic solutions, and other solvents with the diagnostic catheter.

4.3. Preparation

- Apply appropriate anticoagulant therapy to the patient.
- To prevent or reduce clotting and air embolism, always fill the diagnostic catheter with either flushing solution or contrast medium before introducing it into the vascular system. Consider the use of systemic heparinization.
- Use a guidewire with the diagnostic catheter to ensure easy introduction of the catheter.
- The diagnostic catheter can accommodate the maximum guidewire size specified on the label.
- Use coated guidewires with the diagnostic catheter to enhance thromboresistance and lubricity.

4.4. Procedure

- Exercise caution when using 5 Fr and smaller diagnostic catheters and when preoperative imaging indicates that ostial disease is a reasonable expectation. Passing a 5 Fr or smaller diagnostic catheter through an ostial lesion may make detection of ostial disease difficult.
- To avoid damaging the diagnostic catheter after it has been advanced into the vessel, manipulate the guidewire carefully, particularly when negotiating a curve in the diagnostic catheter and when passing through the tip of the diagnostic catheter.
- Reinforce the diagnostic catheter shaft with a guidewire when torque is applied. Excessive torque without the use of a guidewire may result in kinking or twisting, which may make retrieval difficult.
- If resistance is met when advancing or withdrawing the diagnostic catheter, do not proceed until the cause has been determined using fluoroscopic examination.
- Exercise care when removing guidewires from multiple-curve diagnostic catheters.
- Before starting infusion, verify that the catheter has not been kinked or occluded. Failure to do so may cause the catheter to break, rupture, or separate, resulting in damage to the vessel.
- Do not exceed the pressure ratings specified on the label.
- Performance data for the diagnostic catheter is provided in Table 1. Do not exceed the maximum permissible injection pressure.

Instructions for Use English 3

master:20160505, alldocMod:20160504, alldocMod2:20151009, alldocStyle_graph:
20160504, alldocStyle:20160504, master_pageset:20160504, languageLookup:20160330,
table:20160113, characterJoin:20110211, 3of9:20110211, 3of9_svg:20110211

Table 1. Diagnostic Catheter Flow Rate

French Size	Flow rate		Maximum Permissible Pressure (psi/kPa)
	Maximum Flow Rate (ml/sec)	Pressure (psi/kPa)	
5 Fr (1.67 mm)	19	600 psi/4137 kPa	1200 psi/8274 kPa
	27	1200 psi/8274 kPa	
6 Fr (2.00 mm)	28	600 psi/4137 kPa	1200 psi/8274 kPa
	40	1200 psi/8274 kPa	

All measurements were taken using 100 cm length diagnostic catheter samples using a solution of 50% Isovue-300/50% saline by volume at 37°C.

5. Adverse events

Adverse events that may occur or require intervention include, but are not limited to, the following:

- Thrombus formation
- Embolization
- Infection
- Arterial, venous, or heart wall damage or perforation
- Plaque dislodgement
- Hematoma at puncture site
- Cardiac arrhythmia
- Myocardial infarction
- Vascular occlusion
- Death
- Access site pain
- Allergic and neurological reactions
- Vasospasm
- Stroke and Transient Ischemic Attack
- False aneurysm
- Arteriovenous fistula

6. Instructions for use

1. Carefully remove the diagnostic catheter from the package.
2. Flush the diagnostic catheter with heparinized saline solution.
3. The diagnostic catheter may be introduced into the selected vessel by cutdown technique, Seldinger technique, or introducer methods. When using an introducer, ensure that it is a French size greater than or equal to the selected diagnostic catheter.
4. Introduce the diagnostic catheter into the vascular system over a guidewire and carefully advance it under fluoroscopic control.
5. To remove the diagnostic catheter, reinsert the guidewire to straighten the tip of the diagnostic catheter. Slowly withdraw the diagnostic catheter and guidewire under fluoroscopic control. Medtronic recommends using aspiration during withdrawal.

7. Warranty disclaimer

ALTHOUGH THE DIAGNOSTIC CATHETER, HEREFTER REFERRED TO AS "PRODUCT", HAS BEEN MANUFACTURED UNDER CAREFULLY CONTROLLED CONDITIONS, MEDTRONIC, INC., AND AFFILIATES (COLLECTIVELY, "MEDTRONIC") HAVE NO CONTROL OVER CONDITIONS UNDER WHICH THIS PRODUCT IS USED. THE WARNINGS CONTAINED IN THE PRODUCT LABELING PROVIDE MORE DETAILED INFORMATION AND ARE CONSIDERED AN INTEGRAL PART OF THIS DISCLAIMER OF WARRANTY. MEDTRONIC THEREFORE DISCLAIMS ALL WARRANTIES, BOTH EXPRESSED AND IMPLIED, WITH RESPECT TO THE PRODUCT, INCLUDING, BUT NOT LIMITED TO, ANY IMPLIED WARRANTY OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE. SHALL NOT BE LIABLE TO ANY PERSON OR ENTITY FOR ANY MEDICAL EXPENSES OR ANY DIRECT, INCIDENTAL, OR CONSEQUENTIAL DAMAGES CAUSED BY ANY USE, DEFECT, FAILURE, OR MALFUNCTION OF THE PRODUCT, WHETHER A CLAIM FOR SUCH DAMAGES IS BASED UPON WARRANTY, CONTRACT, TORT, OR OTHERWISE. NO PERSON HAS ANY AUTHORITY TO BIND MEDTRONIC TO ANY REPRESENTATION OR WARRANTY WITH RESPECT TO THE PRODUCT.

The exclusions and limitations set out above are not intended to, and should not be construed so as to, contravene mandatory provisions of applicable law. If any part or term of this Disclaimer of Warranty is held to be illegal, unenforceable, or in conflict with applicable law by a court of competent jurisdiction, the validity of the remaining portions of this Disclaimer of Warranty shall not be affected, and all rights and obligations shall be construed and enforced as if this Disclaimer of Warranty did not contain the particular part or term held to be invalid.

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tel 001-763-514-4000
fax 001-763-514-4879

To Whom It May Concern

Medtronic Inc., a corporation organized under the laws of the United States, having its principal place of business at 710 Medtronic Parkway, Minneapolis, MN 55432, USA, as the legal manufacturer of 'DxTerity Diagnostic Angiographic Catheter' hereby confirm that the attached document – «Appendix to the Instruction for use» is a true and original document.

Yours sincerely,

Signed on behalf of Medtronic Inc.
Varsha Kamat
Senior Regulatory Affairs Specialist

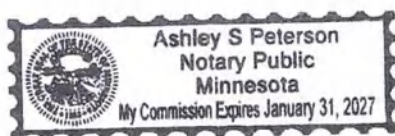
Date: Oct 13, 2022

STATE OF MINNESOTA

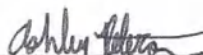
COUNTY OF HENNEPIN

I, Ashley Sophia Peterson, Notary for the State of Minnesota, hereby certify that this document is a true copy of the original filed at Medtronic, Inc.

Dated: October 13, 2022

(stamp)

Signature of Notary Public

Title: Notary Public
My commission expires: January 31, 2027

DxTerity Diagnostic Angiographic Catheter

(When the device is delivered to the Russian Federation, the information contained in this appendix to the technical manual is a priority)

NAME OF MEDICAL DEVICE

DxTerity Diagnostic Angiographic Catheter

1. DxTerity Diagnostic Angiographic Catheter, model: DXT5JL35, DXT6JL35, DXT5JL40, DXT6JL40, DXT5JL40X, DXT6JL40X, DXT5JL45, DXT6JL45, DXT5JL45X, DXT6JL45X, DXT5JL50, DXT6JL50, DXT5JL50X, DXT6JL50X, DXT5JL60, DXT6JL60, DXT5JR35, DXT6JR35, DXT5JR35X, DXT6JR35X, DXT5JR40, DXT6JR40, DXT5JR40X, DXT6JR40X, DXT5JR4ST, DXT6JR4ST, DXT5JR4M, DXT6JR4M, DXT5JR50, DXT6JR50, DXT5JR50X, DXT6JR50X, DXT5JR60, DXT6JR60, DXT5AL10, DXT6AL10, DXT5AL20, DXT6AL20, DXT5AL30, DXT6AL30, DXT5AR10, DXT6AR10, DXT5AR20, DXT6AR20, DXT5ARM, DXT6ARM, DXT5MPA, DXT6MPA, DXT5MPAA, DXT6MPAA, DXT5MPASH, DXT6MPASH, DXT5MPASHA, DXT6MPASHA, DXT5MPB, DXT6MPB, DXT5MPBA, DXT6MPBA, DXT5MPBSH, DXT6MPBSH, DXT5MPBSHA, DXT6MPBSHA, DXT5NOTO, DXT6NOTO, DXT53DRC, DXT63DRC, DXT5RCB, DXT6RCB, DXT5LCB, DXT6LCB, DXT5IMA, DXT6IMA, DXT5PIGSTA, DXT6PIGSTA, DXT5PIGSTX, DXT6PIGSTX, DXT5PIG45A, DXT6PIG45A, DXT5PIG55A, DXT6PIG55A, as part of:

- DxTerity Diagnostic Angiographic Catheter;
- Instructions for Use.

2. DxTerity TRA Diagnostic Angiographic Catheter, model: U5TRAP35, U6TRAP35, U5TRAPSH, U6TRAP35SH, U5TRAP35SHA, U6TRAP35SHA, U5TRAP40, U6TRAP40, U5TRAP40SH, U6TRAP40SH, U5TRAP40SHA, U6TRAP40SHA, U5TRAP45, U6TRAP45, U5TRAP45SH, U6TRAP45SH, U5TRAP45SHA, U6TRAP45SHA, U5TRAN35, U6TRAN35, U5TRAN35SH, U6TRAN35SH, U5TRAN35SHA, U6TRAN35SHA, U5TRAN40, U6TRAN40, U5TRAN40SH, U6TRAN40SH, U5TRAN40SHA, U6TRAN40SHA, U5TRAN45, U6TRAN45, U5TRAN45SH, U6TRAN45SH, U5TRAN45SHA, U6TRAN45SHA, U5TRACK, U6TRACK, U5TRACKSH, U6TRACKSH, U5TRACKSHA, U6TRACKSHA, U5TRACT, U6TRACT, U5TRACTSH, U6TRACTSH, U5TRACTSHA, U6TRACTSHA, as part of:

- DxTerity TRA Diagnostic Angiographic Catheter;
- Instructions for Use.

3. DxTerity Ultra Diagnostic Angiographic Catheter, model: U6ULTRA40SHA, U6ULTRA35SH, U5ULTRA45SHA, U5ULTRA35, U6ULTRA40, U5ULTRA35SH, U6ULTRA45SH, U5ULTRA35SHA, U6ULTRA35, U6ULTRA35SHA, U6ULTRA40SH, U6ULTRA45, U5ULTRA45, U6ULTRA45SHA, U5ULTRA45SH, U5ULTRA40, U5ULTRA40SH, U5ULTRA40SHA, as part of:

- DxTerity Ultra Diagnostic Angiographic Catheter;
- Instructions for Use.

PRESCRIPTION OF MEDICAL DEVICE

It is designed to deliver radiopaque media, guidewires, and therapeutic agents to selected sites in the vascular system.

INFORMATION ABOUT STICKER

The final packaging for the Russian market is accompanied with additional sticker.
Devices supplied in packaging marked with the following possible information:

- Name of medical device in Russian language;
- Name of the manufacturer;
- Country of origin;
- Code UPN (unique product number);
- Code CFN (customer facing number);
- Date and number of Registration Certificate;
- Warning and informing symbols and inscriptions (if applicable)

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TECHNICAL SPECIFICATIONS

Curve Family	Model Number	Product Trade Name	Catheter Length* ²	Catheter French Size	Catheter Curve Type in the package	Device(s) within the package* ¹
Judkins Left	DXT5JL35	DxTerity	100 cm	5	Judkins Left 3.5	5F DxTerity Judkins Left 3.5 Catheter
	DXT6JL35	DxTerity	100 cm	6	Judkins Left 3.5	6F DxTerity Judkins Left 3.5 Catheter
	DXT5JL40	DxTerity	100 cm	5	Judkins Left 4.0	5F DxTerity Judkins Left 4.0 Catheter
	DXT6JL40	DxTerity	100 cm	6	Judkins Left 4.0	6F DxTerity Judkins Left 4.0 Catheter
	DXT5JL40X	DxTerity	125 cm	5	Judkins Left 4.0	5F DxTerity Judkins Left 4.0 Catheter
	DXT6JL40X	DxTerity	125 cm	6	Judkins Left 4.0	6F DxTerity Judkins Left 4.0 Catheter
	DXT5JL45	DxTerity	100 cm	5	Judkins Left 4.5	5F DxTerity Judkins Left 4.5 Catheter

	DXT6JL45	DxTerity	100 cm	6	Judkins Left 4.5	6F DxTerity Judkins Left 4.5 Catheter
	DXT5JL45X	DxTerity	125 cm	5	Judkins Left 4.5	5F DxTerity Judkins Left 4.5 Catheter
	DXT6JL45X	DxTerity	125 cm	6	Judkins Left 4.5	6F DxTerity Judkins Left 4.5 Catheter
	DXT5JL50	DxTerity	100 cm	5	Judkins Left 5.0	5F DxTerity Judkins Left 5.0 Catheter
	DXT6JL50	DxTerity	100 cm	6	Judkins Left 5.0	6F DxTerity Judkins Left 5.0 Catheter
	DXT5JL50X	DxTerity	125 cm	5	Judkins Left 5.0	5F DxTerity Judkins Left 5.0 Catheter
	DXT6JL50X	DxTerity	125 cm	6	Judkins Left 5.0	6F DxTerity Judkins Left 5.0 Catheter
	DXT5JL60	DxTerity	100 cm	5	Judkins Left 6.0	5F DxTerity Judkins Left 6.0 Catheter
	DXT6JL60	DxTerity	100 cm	6	Judkins Left 6.0	6F DxTerity Judkins Left 6.0 Catheter
Judkins Right	DXT5JR35	DxTerity	100 cm	5	Judkins Right 3.5	5F DxTerity Judkins Right 3.5 Catheter
	DXT6JR35	DxTerity	100 cm	6	Judkins Right 3.5	6F DxTerity Judkins Right 3.5 Catheter
	DXT5JR35X	DxTerity	125 cm	5	Judkins Right 3.5	5F DxTerity Judkins Right 3.5 Catheter
	DXT6JR35X	DxTerity	125 cm	6	Judkins Right 3.5	6F DxTerity Judkins Right 3.5 Catheter
	DXT5JR40	DxTerity	100 cm	5	Judkins Right 4.0	5F DxTerity Judkins Right 4.0 Catheter
	DXT6JR40	DxTerity	100 cm	6	Judkins Right 4.0	6F DxTerity Judkins Right 4.0 Catheter
	DXT5JR40X	DxTerity	125 cm	5	Judkins Right 4.0	5F DxTerity Judkins Right 4.0 Catheter
	DXT6JR40X	DxTerity	125 cm	6	Judkins Right 4.0	6F DxTerity Judkins Right 4.0 Catheter

	DXT5JR4ST	DxTerity	100 cm	5	Judkins Right 4.0 ST	5F DxTerity Judkins Right 4.0 ST Catheter
	DXT6JR4ST	DxTerity	100 cm	6	Judkins Right 4.0 ST	6F DxTerity Judkins Right 4.0 ST Catheter
	DXT5JR4M	DxTerity	100 cm	5	Judkins Right 4.0 MOD	5F DxTerity Judkins Right 4.0 MOD Catheter
	DXT6JR4M	DxTerity	100 cm	6	Judkins Right 4.0 MOD	6F DxTerity Judkins Right 4.0 MOD Catheter
	DXT5JR50	DxTerity	100 cm	5	Judkins Right 5.0	5F DxTerity Judkins Right 5.0 Catheter
	DXT6JR50	DxTerity	100 cm	6	Judkins Right 5.0	6F DxTerity Judkins Right 5.0 Catheter
	DXT5JR50X	DxTerity	125 cm	5	Judkins Right 5.0	5F DxTerity Judkins Right 5.0 Catheter
	DXT6JR50X	DxTerity	125 cm	6	Judkins Right 5.0	6F DxTerity Judkins Right 5.0 Catheter
	DXT5JR60	DxTerity	100 cm	5	Judkins Right 6.0	5F DxTerity Judkins Right 6.0 Catheter
	DXT6JR60	DxTerity	100 cm	6	Judkins Right 6.0	6F DxTerity Judkins Right 6.0 Catheter
Amplatz	DXT5AL10	DxTerity	100 cm	5	Amplatz Left 1.0	5F DxTerity Amplatz Left 1.0 Catheter
	DXT6AL10	DxTerity	100 cm	6	Amplatz Left 1.0	6F DxTerity Amplatz Left 1.0 Catheter
	DXT5AL20	DxTerity	100 cm	5	Amplatz Left 2.0	5F DxTerity Amplatz Left 2.0 Catheter
	DXT6AL20	DxTerity	100 cm	6	Amplatz Left 2.0	6F DxTerity Amplatz Left 2.0 Catheter
	DXT5AL30	DxTerity	100 cm	5	Amplatz Left 3.0	5F DxTerity Amplatz Left 3.0 Catheter

	DXT6AL30	DxTernity	100 cm	6	Amplatz Left 3.0	6F DxTernity Amplatz Left 3.0 Catheter
	DXT5AR10	DxTernity	100 cm	5	Amplatz Right 1.0	5F DxTernity Amplatz Right 1.0 Catheter
	DXT6AR10	DxTernity	100 cm	6	Amplatz Right 1.0	6F DxTernity Amplatz Right 1.0 Catheter
	DXT5AR20	DxTernity	100 cm	5	Amplatz Right 2.0	5F DxTernity Amplatz Right 2.0 Catheter
	DXT6AR20	DxTernity	100 cm	6	Amplatz Right 2.0	6F DxTernity Amplatz Right 2.0 Catheter
	DXT5ARM	DxTernity	100 cm	5	Amplatz Right Modified	5F DxTernity Amplatz Right Modified Catheter
	DXT6ARM	DxTernity	100 cm	6	Amplatz Right Modified	6F DxTernity Amplatz Right Modified Catheter
Multipurpose	DXT5MPA	DxTernity	100 cm	5	Multipurpose A	5F DxTernity Multipurpose A Catheter
	DXT6MPA	DxTernity	100 cm	6	Multipurpose A	6F DxTernity Multipurpose A Catheter
	DXT5MPAA	DxTernity	110 cm	5	Multipurpose A	5F DxTernity Multipurpose A Catheter
	DXT6MPAA	DxTernity	110 cm	6	Multipurpose A	6F DxTernity Multipurpose A Catheter
	DXT5MPASH	DxTernity	100 cm	5	Multipurpose A	5F DxTernity Multipurpose A Catheter
	DXT6MPASH	DxTernity	100 cm	6	Multipurpose A	6F DxTernity Multipurpose A Catheter
	DXT5MPASHA	DxTernity	110 cm	5	Multipurpose A	5F DxTernity Multipurpose A Catheter

	DXT6MPASHA	DxTerity	110 cm	6	Multipurpose A	6F DxTerity Multipurpose A Catheter
	DXT5MPB	DxTerity	100 cm	5	Multipurpose B	5F DxTerity Multipurpose B Catheter
	DXT6MPB	DxTerity	100 cm	6	Multipurpose B	6F DxTerity Multipurpose B Catheter
	DXT5MPBA	DxTerity	110 cm	5	Multipurpose B	5F DxTerity™ Multipurpose B Catheter
	DXT6MPBA	DxTerity	110 cm	6	Multipurpose B	6F DxTerity Multipurpose B Catheter
	DXT5MPBSH	DxTerity	100 cm	5	Multipurpose B	5F DxTerity Multipurpose B Catheter
	DXT6MPBSH	DxTerity	100 cm	6	Multipurpose B	6F DxTerity Multipurpose B Catheter
	DXT5MPBSHA	DxTerity	110 cm	5	Multipurpose B	5F DxTerity Multipurpose B Catheter
	DXT6MPBSHA	DxTerity	110 cm	6	Multipurpose B	6F DxTerity Multipurpose B Catheter
Specialty	DXT5NOTO	DxTerity	100 cm	5	No-Torque Right	5F DxTerity No-Torque Right Catheter
	DXT6NOTO	DxTerity	100 cm	6	No-Torque Right	6F DxTerity No-Torque Right Catheter
	DXT53DRC	DxTerity	100 cm	5	3-D Right Coronary	5F DxTerity 3-D Right Coronary Catheter
	DXT63DRC	DxTerity	100 cm	6	3-D Right Coronary	6F DxTerity 3-D Right Coronary Catheter
	DXT5RCB	DxTerity	100 cm	5	Right Coronary Bypass	5F DxTerity Right Coronary Bypass Catheter

	DXT6RCB	DxTernity	100 cm	6	Right Coronary Bypass	6F DxTernity Right Coronary Bypass Catheter
	DXT5LCB	DxTernity	100 cm	5	Left Coronary Bypass	5F DxTernity Left Coronary Bypass Catheter
	DXT6LCB	DxTernity	100 cm	6	Left Coronary Bypass	6F DxTernity Left Coronary Bypass Catheter
	DXT5IMA	DxTernity	100 cm	5	Internal Mammary Artery	5F DxTernity Internal Mammary Artery Catheter
	DXT6IMA	DxTernity	100 cm	6	Internal Mammary Artery	6F DxTernity Internal Mammary Artery Catheter
Pigtails	DXT5PIGSTA	DxTernity	110 cm	5	Pigtail Straight	5F DxTernity Pigtail Straight Catheter
	DXT6PIGSTA	DxTernity	110 cm	6	Pigtail Straight	6F DxTernity Pigtail Straight Catheter
	DXT5PIGSTX	DxTernity	125 cm	5	Pigtail Straight	5F DxTernity Pigtail Straight Catheter
	DXT6PIGSTX	DxTernity	125 cm	6	Pigtail Straight	6F DxTernity Pigtail Straight Catheter
	DXT5PIG45A	DxTernity	110 cm	5	145° Pigtail	5F DxTernity 145° Pigtail Catheter
	DXT6PIG45A	DxTernity	110 cm	6	145° Pigtail	6F DxTernity 145° Pigtail Catheter
	DXT5PIG55A	DxTernity	110 cm	5	155° Pigtail	5F DxTernity 155° Pigtail Catheter

Appendix to the Instruction for use
DxTernity Diagnostic Angiographic Catheter

	DXT6PIG55A	DxTernity	110 cm	6	155° Pigtail	6F DxTernity 155° Pigtail Catheter
Radial	U5TRAP35	DxTernity TRA	100 cm	5	TRApease 3.5	5F DxTernity TRA TRApease 3.5 Catheter
	U6TRAP35	DxTernity TRA	100 cm	6	TRApease 3.5	6F DxTernity TRA TRApease 3.5 Catheter
	U5TRAP35SH	DxTernity TRA	100 cm	5	TRApease 3.5	5F DxTernity TRA TRApease 3.5 Catheter
	U6TRAP35SH	DxTernity TRA	100 cm	6	TRApease 3.5	6F DxTernity TRA TRApease 3.5 Catheter
	U5TRAP35SHA	DxTernity TRA	110 cm	5	TRApease 3.5	5F DxTernity TRA TRApease 3.5 Catheter
	U6TRAP35SHA	DxTernity TRA	110 cm	6	TRApease 3.5	6F DxTernity TRA TRApease 3.5 Catheter
	U5TRAP40	DxTernity TRA	100 cm	5	TRApease 4.0	5F DxTernity TRA TRApease 4.0 Catheter
	U6TRAP40	DxTernity TRA	100 cm	6	TRApease 4.0	6F DxTernity TRA TRApease 4.0 Catheter
	U5TRAP40SH	DxTernity TRA	100 cm	5	TRApease 4.0	5F DxTernity TRA TRApease 4.0 Catheter
	U6TRAP40SH	DxTernity TRA	100 cm	6	TRApease 4.0	6F DxTernity TRA TRApease 4.0 Catheter
	U5TRAP40SHA	DxTernity TRA	110 cm	5	TRApease 4.0	5F DxTernity TRA TRApease 4.0 Catheter
	U6TRAP40SHA	DxTernity TRA	110 cm	6	TRApease 4.0	6F DxTernity TRA TRApease 4.0 Catheter

	U5TRAP45	DxTerity TRA	100 cm	5	TRApease 4.5	5F DxTerity TRA TRApease 4.5 Catheter
	U6TRAP45	DxTerity TRA	100 cm	6	TRApease 4.5	6F DxTerity TRA TRApease 4.5 Catheter
	U5TRAP45SH	DxTerity TRA	100 cm	5	TRApease 4.5	5F DxTerity TRA TRApease 4.5 Catheter
	U6TRAP45SH	DxTerity TRA	100 cm	6	TRApease 4.5	6F DxTerity TRA TRApease 4.5 Catheter
	U5TRAP45SHA	DxTerity TRA	110 cm	5	TRApease 4.5	5F DxTerity TRA TRApease 4.5 Catheter
	U6TRAP45SHA	DxTerity TRA	110 cm	6	TRApease 4.5	6F DxTerity TRA TRApease 4.5 Catheter
	U5TRAN35	DxTerity TRA	100 cm	5	TRAnsformer 3.5	5F DxTerity TRA TRAnsformer 3.5 Catheter
	U6TRAN35	DxTerity TRA	100 cm	6	TRAnsformer 3.5	6F DxTerity TRA TRAnsformer 3.5 Catheter
	U5TRAN35SH	DxTerity TRA	100 cm	5	TRAnsformer 3.5	5F DxTerity TRA TRAnsformer 3.5 Catheter
	U6TRAN35SH	DxTerity TRA	100 cm	6	TRAnsformer 3.5	6F DxTerity TRA TRAnsformer 3.5 Catheter
	U5TRAN35SHA	DxTerity TRA	110 cm	5	TRAnsformer 3.5	5F DxTerity TRA TRAnsformer 3.5 Catheter
	U6TRAN35SHA	DxTerity TRA	110 cm	6	TRAnsformer 3.5	6F DxTerity TRA TRAnsformer 3.5 Catheter
	U5TRAN40	DxTerity TRA	100 cm	5	TRAnsformer 4.0	5F DxTerity TRA

						TRAnsformer 4.0 Catheter
	U6TRAN40	DxTerity TRA	100 cm	6	TRAnsformer 4.0	6F DxTerity TRA TRAnsformer 4.0 Catheter
	U5TRAN40SH	DxTerity TRA	100 cm	5	TRAnsformer 4.0	5F DxTerity TRA TRAnsformer 4.0 Catheter
	U6TRAN40SH	DxTerity TRA	100 cm	6	TRAnsformer 4.0	6F DxTerity TRA TRAnsformer 4.0 Catheter
	U5TRAN40SHA	DxTerity TRA	110 cm	5	TRAnsformer 4.0	5F DxTerity TRA TRAnsformer 4.0 Catheter
	U6TRAN40SHA	DxTerity TRA	110 cm	6	TRAnsformer 4.0	6F DxTerity TRA TRAnsformer 4.0 Catheter
	U5TRAN45	DxTerity TRA	100 cm	5	TRAnsformer 4.5	5F DxTerity TRA TRAnsformer 4.5 Catheter
	U6TRAN45	DxTerity TRA	100 cm	6	TRAnsformer 4.5	6F DxTerity TRA TRAnsformer 4.5 Catheter
	U5TRAN45SH	DxTerity TRA	100 cm	5	TRAnsformer 4.5	5F DxTerity TRA TRAnsformer 4.5 Catheter
	U6TRAN45SH	DxTerity TRA	100 cm	6	TRAnsformer 4.5	6F DxTerity TRA TRAnsformer 4.5 Catheter
	U5TRAN45SHA	DxTerity TRA	110 cm	5	TRAnsformer 4.5	5F DxTerity TRA TRAnsformer 4.5 Catheter
	U6TRAN45SHA	DxTerity TRA	110 cm	6	TRAnsformer 4.5	6F DxTerity TRA TRAnsformer 4.5 Catheter
	U5TRACK	DxTerity TRA	100 cm	5	TRAckeR	5F DxTerity TRA TRAckeR Catheter

	U6TRACK	DxTernity TRA	100 cm	6	TRAcKer	6F DxTernity TRA TRAcKer Catheter
	U5TRACKSH	DxTernity TRA	100 cm	5	TRAcKer	5F DxTernity TRA TRAcKer Catheter
	U6TRACKSH	DxTernity TRA	100 cm	6	TRAcKer	6F DxTernity TRA TRAcKer Catheter
	U5TRACKSHA	DxTernity TRA	110 cm	5	TRAcKer	5F DxTernity TRA TRAcKer Catheter
	U6TRACKSHA	DxTernity TRA	110 cm	6	TRAcKer	6F DxTernity TRA TRAcKer Catheter
	U5TRACT	DxTernity TRA	100 cm	5	TRActive	5F DxTernity TRA TRActive Catheter
	U6TRACT	DxTernity TRA	100 cm	6	TRActive	6F DxTernity TRA TRActive Catheter
	U5TRACTSH	DxTernity TRA	100 cm	5	TRActive	5F DxTernity TRA TRActive Catheter
	U6TRACTSH	DxTernity TRA	100 cm	6	TRActive	6F DxTernity TRA TRActive Catheter
	U5TRACTSHA	DxTernity TRA	110 cm	5	TRActive	5F DxTernity TRA TRActive Catheter
	U6TRACTSHA	DxTernity TRA	110 cm	6	TRActive	6F DxTernity TRA TRActive Catheter
Ultra	U6ULTRA40SHA	DxTernity TRA	110 cm	6	ULTRA	6F DxTernity TRA
	U6ULTRA35SH	DxTernity TRA	100 cm	6	ULTRA	6F DxTernity TRA
	U5ULTRA45SHA	DxTernity TRA	110 cm	5	ULTRA	5F DxTernity TRA

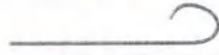
Appendix to the Instruction for use
DxTerity Diagnostic Angiographic Catheter

	U5ULTRA35	DxTerity TRA	100 cm	5	ULTRA	5F DxTerity TRA
	U6ULTRA40	DxTerity TRA	100 cm	6	ULTRA	6F DxTerity TRA
	U5ULTRA35SH	DxTerity TRA	100 cm	5	ULTRA	5F DxTerity TRA
	U6ULTRA45SH	DxTerity TRA	100 cm	6	ULTRA	6F DxTerity TRA
	U5ULTRA35SHA	DxTerity TRA	110 cm	5	ULTRA	5F DxTerity TRA
	U6ULTRA35	DxTerity TRA	100 cm	6	ULTRA	6F DxTerity TRA
	U6ULTRA35SHA	DxTerity TRA	110 cm	6	ULTRA	6F DxTerity TRA
	U6ULTRA40SH	DxTerity TRA	100 cm	6	ULTRA	6F DxTerity TRA
	U6ULTRA45	DxTerity TRA	100 cm	6	ULTRA	6F DxTerity TRA
	U5ULTRA45	DxTerity TRA	100 cm	5	ULTRA	5F DxTerity TRA
	U6ULTRA45SHA	DxTerity TRA	110 cm	6	ULTRA	6F DxTerity TRA
	U5ULTRA45SH	DxTerity TRA	100 cm	5	ULTRA	5F DxTerity TRA
	U5ULTRA40	DxTerity TRA	100 cm	5	ULTRA	5F DxTerity TRA
	U5ULTRA40SH	DxTerity TRA	100 cm	5	ULTRA	5F DxTerity TRA
	U5ULTRA40SHA	DxTerity TRA	110 cm	5	ULTRA	5F DxTerity TRA

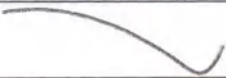
*¹: 5F(5)= 1,67 mm; 6F(6)= 2 mm; permissible deviation $\pm 5\%$

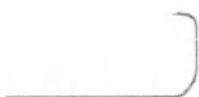
*²: permissible deviation ± 2 cm.

List of distal curve shapes - DxTerity Catheters

Curve Family	Curve Description	Curve Picture	Models
Judkins Left	Judkins Left 3.5		DXT5JL35 DXT6JL35
	Judkins Left 4.0		DXT5JL40 DXT6JL40 DXT5JL40X DXT6JL40X
	Judkins Left 4.5		DXT5JL45 DXT6JL45 DXT5JL45X DXT6JL45X
	Judkins Left 5.0		DXT5JL50 DXT6JL50 DXT5JL50X DXT6JL50X
	Judkins Left 6.0		DXT5JL60 DXT6JL60

Curve Family	Curve Description	Curve Picture	Models
Judkins Right	Judkins Right 3.5		DXT5JR35 DXT6JR35 DXT5JR35X DXT6JR35X
	Judkins Right 4.0		DXT5JR40 DXT6JR40 DXT5JR40X DXT6JR40X
	Judkins Right 4.0ST		DXT5JR4ST DXT6JR4ST
	Judkins Right 4.0M		DXT5JR4M DXT6JR4M
	Judkins Right 5.0		DXT5JR50 DXT6JR50 DXT5JR50X DXT6JR50X
	Judkins Right 6.0		DXT5JR60 DXT6JR60
Amplatz	Amplatz Left 1.0		DXT5AL10 DXT6AL10
	Amplatz Left 2.0		DXT5AL20 DXT6AL20
	Amplatz Left 3.0		DXT5AL30 DXT6AL30
	Amplatz Right 1.0		DXT5AR10 DXT6AR10
	Amplatz Right 2.0		DXT5AR20 DXT6AR20
	Amplatz Right Modified		DXT5ARM DXT6ARM
Multipurpose	Multipurpose A		DXT5MPA DXT6MPA DXT5MPAA DXT6MPAA DXT5MPASH DXT6MPASH DXT5MPASHA DXT6MPASHA
	Multipurpose B		DXT5MPB DXT6MPB DXT5MPBA DXT6MPBA DXT5MPBSH DXT6MPBSH DXT5MPBSHA DXT6MPBSHA
Specialty	No-Torque Right		DXT5NOTO DXT6NOTO

Curve Family	Curve Description	Curve Picture	Models
	3-D Right Coronary		DXT53DRC DXT63DRC
	Right Coronary Bypass		DXT5RCB DXT6RCB
	Left Coronary Bypass		DXT5LCB DXT6LCB
	Internal Mammary Artery		DXT5IMA DXT6IMA
Pigtail	Pigtail Straight		DXT5PIGSTA DXT6PIGSTA DXT5PIGSTX DXT6PIGSTX
	145° Pigtail		DXT5PIG45A DXT6PIG45A
	155° Pigtail		DXT5PIG55A DXT6PIG55A
Radial Curves	TRapease 3.5		U5TRAP35 U6TRAP35 U5TRAP35SH U6TRAP35SH U5TRAP35SHA U6TRAP35SHA
	TRapease 4.0		U5TRAP40 U6TRAP40 U5TRAP40SH U6TRAP40SH U5TRAP40SHA U6TRAP40SHA
	TRapease 4.5		U5TRAP45 U6TRAP45 U5TRAP45SH U6TRAP45SH U5TRAP45SHA U6TRAP45SHA
	TRAnsformer 3.5		U5TRAN35 U6TRAN35 U5TRAN35SH U6TRAN35SH U5TRAN35SHA U6TRAN35SHA
	TRAnsformer 4.0		U5TRAN40 U6TRAN40 U5TRAN40SH U6TRAN40SH U5TRAN40SHA U6TRAN40SHA
	TRAnsformer 4.5		U5TRAN45 U6TRAN45 U5TRAN45SH U6TRAN45SH

Curve Family	Curve Description	Curve Picture	Models
			U5TRAN45SHA U6TRAN45SHA
	TRacker		U5TRACK U6TRACK U5TRACKSH U6TRACKSH U5TRACKSHA U6TRACKSHA
	TRActive		U5TRACT U6TRACT U5TRACTSH U6TRACTSH U5TRACTSHA U6TRACTSHA
Ultra Curves	Ultra 3.5		U6ULTRA35SH U5ULTRA35 U5ULTRA35SH U5ULTRA35SHA U6ULTRA35 U6ULTRA35SHA
	Ultra 4.0		U6ULTRA40SHA U6ULTRA40 U6ULTRA40SH U5ULTRA40 U5ULTRA40SH U5ULTRA40SHA
	Ultra 4.5		U5ULTRA45SHA U6ULTRA45SH U6ULTRA45 U5ULTRA45 U6ULTRA45SHA U5ULTRA45SH

Coronary catheters include Judkins Left, Judkins Right, Amplatz, Multipurpose, Specialty, Radial Curves, Ultra families

Coronary Catheter

PRODUCT CHARACTERISTIC	SPECIFICATION**
Effective Length	100 cm ± 2 cm
	110 cm ± 2 cm
	125 cm ± 2 cm
Weight	16 g ± 2 g
Soft Tip ID	5F = 0.042" (1.07 mm) minimum
	6F = 0.052" (1.32 mm) minimum
Shaft OD	5F = 0.070" (1.78 mm) maximum
	6F = 0.081" (2.06 mm) maximum
Soft Tip OD	5F = 0.072" (1.83 mm) maximum
	6F = 0.083" (2.11 mm) maximum
Soft Tip Tensile	1.0 lbf (4.5N) minimum
Width of the blue polymer layer	4 mm ± 0.4 mm
Hub/Shaft Tensile	5F = 2.2 lbf (10N) minimum
	6F = 3.4 lbf (15N) minimum

PRODUCT CHARACTERISTIC	SPECIFICATION**
Hub/Shaft Leak/Burst	No leaks/burst up to 1200 psi Minimum
Shaft Stiffness	5F: 17% $\pm$ 7% of a full-scale bending moment of 0.35 in-lbf at a deflection angle of 40° 6F: 30% $\pm$ 10% of a full-scale bending moment of 0.35 in-lbf at a deflection angle of 40°
Shaft Kink	60° minimum
Straight Line Torque	5F = 0.8 oz-in minimum at 40° angular rotation 6F = 1.2 oz-in minimum at 40° angular rotation
Torque to Kink	5F = 245° minimum 6F = 175° minimum
Curve Retention (at body temperature)	Primary: 128° maximum Secondary: 68° maximum
Primary Curve Draw Through Stiffness (at body temperature)	5F = 0.16 lbf +/- 0.10 lbf 6F = 0.26 lbf +/- 0.14 lbf
Secondary Curve Draw Through Stiffness (at body temperature)	5F = 0.16 lbf +/- 0.10 lbf 6F = 0.26 lbf +/- 0.14 lbf
Soft Tip Compression	5F = 40 gf +/- 25 gf 6F = 60 gf +/- 35 gf
Surface Lubricity (at body temperature)	Lubricity Maximum 500 gf (Avg runs 1-3)
Flow Rate	Flow Rate at 600 psi Max 5F: 15 ml/sec. minimum / 19 ml/sec. maximum 6F: 17 ml/sec. minimum / 28 ml/sec. maximum Flow Rate at 1200 psi Max 5F: 21 ml/sec. minimum / 27 ml/sec. maximum 6F: 25 ml/sec. minimum / 40 ml/sec. maximum 50/50 Isovue-300 contrast/saline mixture @ 37°C
Luer Taper	Luer taper must meet ISO 594-1/ISO 594-2 (ISO 80369-7)
Cannula size, L x W x T, mm $\pm$ 5%	19 x 15 x 5
Tightness	No air leakage at the connection nodes and the connection of the outlet with Luer at a load of 700 kPa
Radiographic	Visibility Medical device is radiopaque
** 5F = 1.67 mm; 6 F= 2 mm.	

Pigtail Catheter Physical

PRODUCT CHARACTERISTIC	SPECIFICATION
Effective Length	110 cm $\pm$ 2cm 125 cm $\pm$ 2cm
Weight	16 g $\pm$ 2 g
Tip ID	0.040" (1.02mm) minimum
Shaft OD	Pigtail Catheter utilizes the same braided shaft as Coronary Catheter. See Coronary Catheter Shaft OD specification
Segment Joint OD	5F= 0.072" (1.83mm) maximum 6F= 0.083" (2.11mm) maximum
Segment Joint Tensile	5F = 2.2 lbf (10N) minimum 6F = 3.4 lbf (15N) minimum
Segment Joint Flex	No complete material fractures through the wall at or below 40 flexures for each sample

PRODUCT CHARACTERISTIC	SPECIFICATION
Hub/Shaft Tensile Strength	Pigtail Catheter utilizes the same Hub and braided Shaft as Coronary Catheter. See Coronary Catheter Hub/Shaft Tensile Strength specification
Hub/Shaft Leak/Burst	No leaks/burst up to 1200 psi Minimum
Shaft Stiffness	Pigtail Catheter utilizes the same braided Shaft as Coronary Catheter. See Coronary Catheter Shaft Stiffness specification
Shaft Kink	Pigtail Catheter utilizes the same braided Shaft as Coronary Catheter. See Coronary Catheter Shaft Kink specification
Straight Line Torque	Pigtail Catheter utilizes the same braided Shaft as Coronary Catheter. See Coronary Catheter Straight Line Torque specification
Torque to Kink	Pigtail Catheter utilizes the same braided Shaft as Coronary Catheter. See Coronary Catheter Torque to Kink specification
Pigtail Curve Draw Through Stiffness (at body temperature)	0.22 lbf +/- 0.15 lbf
Surface Lubricity	Pigtail Catheter utilizes the same braided shaft as Coronary Catheter. See Coronary Catheter Surface Lubricity specification
Flow Rate	Flow Rate at 600 psi Max 5F: 15 ml/sec. minimum 6F: 17 ml/sec. minimum Flow Rate at 1200 psi Max 5F: 21 ml/sec. minimum 6F: 25 ml/sec. minimum 50/50 Isovue-300 contrast/saline mixture @ 37C
Luer Taper	Pigtail Catheter utilizes the same Hub as Coronary Catheter. See Coronary Catheter Luer Taper specification
Cannula size, l x W x T, mm $\pm$ 5%	19 x 15 x 5
Number of holes at the distal end	6
Hole diameter, mm, $\pm$ 5 %	0.9
Tightness	No air leakage at the connection nodes and the connection of the outlet with Luer at a load of 700 kPa
Radiographic	Visibility Medical device is radiopaque
** 5F = 1,67 mm; 6 F= 2 mm.	

STERILIZATION

This device is single-use.

Device is sterile and non-pyrogenic.

Sterilized by EtO.

Sterility Assurance Level (SAL) - 10^{-6}

Term preservation of sterility during 3 years (for Coronary).

Term preservation of sterility during 2 years (for Pigtail).

Primary packaging in accordance with ISO 11607.

CONDITIONS OF TRANSPORTATION AND STORAGE

Store diagnostic catheters in a cool, dry, and dark area.

Do not expose diagnostic catheters to UV light.

Do not heat or bend the diagnostic catheter before use.

Temperature from + 20°C to +25°C

Relative Humidity of 20% -80%

Pressure from 100 kPa to 103 kPa

OPERATING CONDITIONS

Do not exceed the pressure ratings specified on the label.

Temperature, °C from +32 to + 42

Humidity, % Up To 100

Relative pressure, kPa from 700 to 1060

SHELF LIFE

Coronary - 3 years.

Pigtail - 2 years.

PROCEDURE AND CONDITIONS FOR THE DISPOSAL OR DESTRUCTION OF A MEDICAL DEVICE

Devices should be utilized according to SanPiN 2.1.3684-21 "Sanitary and epidemiological requirements for the maintenance of territories of urban and rural settlements, for water bodies, drinking water and drinking water supply, atmospheric air, soils, residential premises, operation of industrial, public premises, organization and implementation of sanitary and anti-epidemic (preventive) measures)".

After use, this product may be a potential biohazard (hazard class B). Handle and dispose of all such device in accordance with accepted medical practice and applicable local, regional and national laws and regulations.

THE LIST OF NATIONAL AND INTERNATIONAL REGULATORY DOCUMENTS/STANDARDS, WHICH MATCHES MEDICAL DEVICE

EN ISO 13485, ISO 10555-1, ISO 594-2, ISO 14644-1, ISO 14644-2, MDD 93/42/EEC, EN ISO 10993-1, EN ISO 10993-4, EN ISO 10993-5, EN ISO 10993-10, EN ISO 10993-11, EN ISO 10993-12, EN 556-1, EN ISO 10993-7, ISO 11135, ANSI/AAMI ST72, EN ISO 11607-1, EN ISO 11607-2, EN 980/ISO 15223-1, EN 1041, MEDDEV.2.7.1, EN ISO 14971.

CUSTOMER SERVICE

Contact for all questions, including complaints of medical devices the Authorized Representative in Russia:

Medtronic, LLC,

123112, Russia, Moscow, in.ter.c. Presnensky Municipal District, Presnenskaya nab., 10, floor 9, office III, room 41.

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© 2016 Medtronic