Medtronic**ARCTIC FRONT ADVANCE™ 2AF233, 2AF283**

Cardiac Cryoablation Catheter
 Srdečni katétr pro kryoablaci
 Kardielt kryoablationskateter
 Kardialer Kryoablationskatheter
 Καρδιακός καθετήρας κρυοκατάλυσης
 Catéter de crioablación cardíaca
 Cathéter de cryoablation cardiaque
 Kardiális krioablációs katéter
 Catetere per crioablazione cardiaca
 Cardiale cryoablatiekatheter
 Kardialt kryoablasijskateter
 Cewnik do krioablacji serca
 Cateter cardiaco de crioablación
 Катетер для внутрисердечной кrioабляции
 Kardiálny kryoablačný katéter
 Cardiac Cryoablation-kateter
 Kardiyak Cryoablasyon Kateteri

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True and exact copy of the original document on file at
 Medtronic, Inc.

Robyn VanReese, Principal Regulatory Affairs Specialist
 August 21, 2020

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CE
 0123

The following list includes trademarks or registered trademarks of Medtronic in the United States and possibly in other countries. All other trademarks are the property of their respective owners.

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Achieve, Arctic Front, Arctic Front Advance, CryoConsole, FlexCath, Medtronic

Explanation of symbols • Vysvětlení symbolů • Symbolforklaring • Erläuterung der Symbole • Επεξήγηση των συμβόλων • Explicación de los símbolos • Explication des symboles • A szimbólumok magyarázata • Spiegazione dei simboli • Verklaring van de symbolen • Symbolforklaring • Objasnienie symboli • Explicação dos símbolos • Объяснение символов • Vysvetlivky k symbolom • Förklaring av symboler • Sembollerin açıklaması

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1 Description


The Arctic Front Advance Cardiac Cryoablation Catheter (the catheter or the Arctic Front Advance Cryoballoon) is a flexible, over-the-wire balloon catheter used to ablate cardiac tissue. It is used together with a compatible Medtronic 12 Fr inner diameter sheath (the sheath), the CryoConsole, and related components. For device compatibility questions, contact Medtronic Technical Support.

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Figure 1. Arctic Front Advance Cardiac Cryoablation Catheter

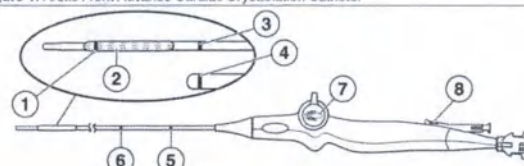- | | |
|---------------------------------------|----------------------------------|
| 1 Distal catheter radiopaque marker | 5 Proximal catheter shaft marker |
| 2 Deflated balloon segment | 6 Distal catheter shaft marker |
| 3 Proximal catheter radiopaque marker | 7 Deflection mechanism |
| 4 FlexCath sheath radiopaque marker | 8 Blue push button |

The Arctic Front Advance Cardiac Cryoablation Catheter is available in 2 models, as described in the following table:

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- 1 Arctic Front Advance Cardiac Cryoablation Catheter
- product documentation

2 Indications for use

The Arctic Front Advance Cardiac Cryoablation Catheter is indicated for the treatment of patients with atrial fibrillation.

3 Contraindications

The Arctic Front Advance Cardiac Cryoablation Catheter is contraindicated as follows:

- in the ventricle because of the danger of catheter entrapment in the chordae tendineae
- in patients with active systemic infections
- in conditions where the manipulation of the catheter within the heart would be unsafe (for example, intracardiac mural thrombus)
- in patients with cryoglobulinemia
- in patients with one or more pulmonary vein stents

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- Do not inflate the balloon while the catheter is positioned inside a pulmonary vein. Always inflate the balloon in the atrium and then position it at the pulmonary vein ostium. Inflating the balloon in the pulmonary vein may result in vascular injury.
- If the balloon cannot be inflated or deflated using the CryoConsole, have a Manual Retraction Kit on hand during the procedure. (Refer to the *CryoConsole Operator's Manual* for more detailed instructions on the Manual Retraction Kit).

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- Do not use excessive force to advance, withdraw, or rotate the catheter, especially if resistance is encountered. Excessive force may lead to catheter damage, including kinking of the guide wire lumen within the balloon segment.
- Do not use the catheter if it is kinked, damaged, or cannot be straightened.
- Straighten the shaft before inserting or withdrawing the catheter.
- Do not at any time preshape or bend the catheter shaft or balloon segment. Bending or kinking the catheter shaft may damage internal structures and increase the risk of catheter failure. Prebending of the distal curve may damage the catheter.
- Catheter advancement should be performed using fluoroscopy or other appropriate techniques.
- Do not position the cryoballoon catheter within the tubular portion of the pulmonary vein to minimize phrenic nerve injury and pulmonary vein stenosis.

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Do not resterilize – Do not resterilize this device for the purpose of reuse. Resterilization may compromise the structural integrity of the device or create a risk of contamination from the device that could result in patient injury, illness, or death.

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Environmental limits – Perform cryoablation procedures only within the environmental parameters. Operating outside these parameters may prevent the start or completion of a cryoablation procedure. Refer to Chapter 7, "Specifications", page 8 for environmental parameters.

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Guide wire compatibility – Use only 0.081 cm (0.032 in) or 0.089 cm (0.035 in) guide wires with the catheter. Using another guide wire may damage the catheter or compromise the procedure.

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Pulmonary vein narrowing or stenosis – Catheter ablation procedures inside or near pulmonary veins may induce pulmonary vein narrowing or stenosis. Do not ablate in the tubular portion of the pulmonary vein. The occurrence of this complication may necessitate percutaneous angioplasty or surgical intervention.

Qualified users – This cryoablation system should be used only by or under the supervision of physicians trained in cryoablation procedures.

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Sterile package inspection – Inspect the sterile packaging and catheter before use. If the sterile packaging or catheter is damaged, do not use the catheter. Contact your Medtronic representative.

System compatibility – Use only Medtronic cryoablation catheters, refrigerant tanks, and components with the CryoConsole. The safety and use of other catheters or components has not been tested.

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- Access site complications (e.g. bruising, ecchymosis)
- Anemia
- Anxiety
- Arrhythmia (e.g. atrial flutter, bradycardia, heart block, tachycardia)
- Back pain
- Bleeding from puncture sites
- Bronchial constriction
- Bronchial fistula
- Bronchitis
- Bruising
- Cardiac tamponade
- Cardiopulmonary arrest
- Cerebral vascular accident
- Chest discomfort/pain/pressure
- Cold feeling
- Coronary artery spasm
- Cough
- Death
- Diarrhea
- Dizziness
- Embolism
- Esophageal damage (including atrio-esophageal fistula)
- Fatigue
- Fever
- Headache
- Hemoptysis
- Hypotension/hypertension
- Infection (e.g. pericarditis, sepsis, urinary)
- Lightheadedness
- Myocardial infarction
- Nausea/vomiting
- Perforation
- Pericardial effusion
- Phrenic nerve injury
- Pleural effusion
- Pneumonia
- Pneumothorax
- Pseudoaneurysm
- Pulmonary edema
- Pulmonary hemorrhage
- Pulmonary vein dissection
- Pulmonary vein stenosis
- Shivering
- Shortness of breath
- Sore throat
- Transient ischemic attack
- Vagal nerve injury (e.g. gastroparesis)
- Vasovagal reaction
- Visual changes (e.g. blurred vision)

6 Instructions for use

6.1 Connecting the Arctic Front Advance Cardiac Cryoablation Catheter

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1. Connect the auto connection box to the CryoConsole.
Note: The ECG cable is not required for an Arctic Front Advance Cryoballoon procedure, and should not be connected to the auto connection box.
2. Connect the catheter to a sterile electrical umbilical cable and then pass the other end of the cable out of the sterile field and connect it to the auto connection box.
3. Connect the catheter to a sterile coaxial umbilical cable and then pass the other end of the cable out of the sterile field and connect it to the CryoConsole.
4. Enable vacuum on CryoConsole.

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Notes:

- Before introducing the catheter into the patient, test the deflection mechanism on the handle to ensure it is operational.
 - Always use the deflection mechanism on the handle to straighten the shaft before insertion or withdrawal of the catheter.
 - The Achieve family of mapping catheters is compatible for use with the Medtronic Arctic Front Advance Cryoablation Catheter and may be used to support and position the catheter.
1. Using aseptic technique, create a vascular access with an appropriate introducer. Obtain left atrial transeptal access using a transeptal sheath, its dilator, and needle.
 - Place standard diagnostic pacing catheters.
 - Visualize left atrial anatomy to help select a balloon size. Select balloon size.
 - Selection of balloon size should be based on pulmonary vein diameter and shape, the surrounding anatomy, and desired position of the balloon outside the tubular portion of the pulmonary vein.
 2. Remove the transeptal sheath and dilator, leaving the guide wire positioned preferably in the left superior pulmonary vein.
 3. Advance the sheath and dilator over the wire into the left atrium.
 4. Slowly remove the guide wire and dilator from the sheath.
 5. Aspirate and flush the sheath.
 6. Obtain a catheter. Connect the Tuohy-Borst adapter and manifold to the push button Luer.
 7. Flush the guide wire lumen, ensuring there are no air bubbles in the guide wire lumen or connections.
 8. Load the circular mapping catheter or guide wire into the guide wire lumen.
 9. Pull the circular mapping catheter or guide wire back into the balloon catheter until it is just inside the balloon catheter tip.
 10. Flush the guide wire lumen again, ensuring there are no air bubbles in the guide wire lumen or connections.
 11. Retract the sleeve from the balloon onto the catheter shaft.
 12. Submerge and rinse the balloon to remove trapped air bubbles.
 13. Advance the sleeve back over the balloon while the balloon is still submerged.
 14. Place the sleeve adjacent to the hemostasis valve of the sheath.
Note: To avoid damaging the valve or introducing air, do not push the sleeve through the valve opening.
 15. Insert the catheter into the sheath.
Note: Follow institutional procedures to aspirate and flush the sheath.
 16. Advance the balloon until the distal tip of the catheter aligns with the distal tip of the sheath using fluoroscopic guidance or other appropriate visualization techniques.
 17. Advance the circular mapping catheter or guide wire to the target pulmonary vein using fluoroscopic guidance or other appropriate visualization techniques.
 18. Advance the balloon over the circular mapping catheter or guide wire into the left atrium using fluoroscopic guidance or other appropriate visualization techniques.
 19. Inflate the balloon in the left atrium by pressing the Start button on the CryoConsole control panel for 2 seconds.

20. Position the catheter at the ostium of the target pulmonary vein and not inside the tubular portion of the pulmonary vein.
21. Verify the balloon position by injecting a mixture of 50/50 contrast/saline into the catheter guide wire lumen port or by using other appropriate techniques. Be sure to flush the guide wire lumen with saline after each contrast injection.
- Notes:**
- To improve balloon position and support, reposition the circular mapping catheter. If a stable balloon position cannot be obtained, exchange the circular mapping catheter for a guide wire.
 - Before exchanging the circular mapping catheter for a guide wire, retract the catheter into the sheath.
 - Follow contrast labeling and institutional procedures regarding the appropriate medical strategies to minimize the risk to the patient associated with using contrast.
22. Perform the cryoablation.
- Notes:**
- During the initial ablation phase, monitor the balloon's position using appropriate visualization techniques. Moisture in the system may cause the inflated balloon to deflate and then re-inflate during the initial transition phases.
 - Set the ablation duration on the CryoConsole screen.
 - Physicians may modify the preset ablation duration based on clinical judgment. It is not recommended that the ablation duration exceeds 240 seconds.
 - The balloon's outer diameter varies from inflation to cryoablation and may cause the balloon to shift.
23. Wait for the cryoablation phase to complete (at the end of the preset duration). The balloon remains inflated and the thawing phase begins.
- Note:** At any time, the ablation can be stopped by pressing the Stop Current Action button on the CryoConsole control panel or by pressing the Stop CryoAblation button on the screen.
24. During the thawing phase, observe the temperature indicator on the screen. The balloon deflates automatically when the temperature reaches 20°C.
25. Before moving the balloon, use appropriate techniques to ensure that the catheter is not adhered to tissue.
26. Determine effective ablation of the cardiac tissue by assessing electrical isolation of the pulmonary vein from the left atrium (entrance and exit block) after the cryoablation is complete.
- Note:** As needed, perform additional treatments by positioning the balloon differently in the same pulmonary vein.
27. Position the catheter at the ostium of the next target pulmonary vein using the sheath, circular mapping catheter, or guide wire. Return to Step 17 and continue ablation.
28. To retract the balloon into the sheath, perform the following steps:
- Note:** Ensure that the distal tip of the catheter is free to move to its maximum length.
- Use the deflection mechanism on the sheath handle to straighten the sheath. Use the deflection mechanism on the catheter handle to straighten the catheter.
 - Inflate the balloon.
 - Perform the following two steps simultaneously:
Advance the blue push button on the catheter handle, as shown in Figure 2. This causes the balloon to extend to maximum length and wrap tightly.
Deflate the balloon by pressing the Stop Current Action button on the CryoConsole control panel, or by selecting the Deflate Balloon option on the screen.
 - Retract the catheter into the sheath.

Figure 2. Catheter handle

29. Remove the catheter from the patient.
30. After the procedure, follow the instructions to "Shut down the system" in the *CryoConsole Operator's Manual*. Follow the prompts to close the refrigerant tank and replace the cap on the coaxial port. Confirm these actions when prompted by the CryoConsole.

7 Specifications

Catheter shaft size	3.5 mm (10.5 Fr; 0.14 in)
Tip outer diameter	3 mm (9 Fr; 0.12 in)
Recommended introducer sheath	compatible Medtronic 12 Fr inner diameter sheath
Inner diameter of guide wire lumen	1.27 mm (0.05 in) nominal
Inflated balloon diameter	2AF233 – 23 mm (0.91 in) 2AF283 – 28 mm (1.10 in)
Effective length (with balloon inflated)	95.0 ±2.0 cm (37.40 ±0.80 in)
Number of thermocouples	1
Environmental parameters	
Storage temperature	15°C to 30°C (59°F to 86°F)
Transit temperature	–35°C to 45°C (–31°F to 113°F); up to 85% relative humidity (non-condensing)
Operation	15°C to 30°C (59°F to 86°F) at altitudes less than 2400 m (8000 feet) above sea level

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English

1 Description


The Arctic Front Advance cardiac cryoablation catheter (the catheter or the Arctic Front Advance Cryoballoon) is a flexible, over-the-wire balloon catheter used to ablate cardiac tissue. The catheter delivers cryoenergy at the antrum of the pulmonary veins to create circumferential lesions in order to ablate cardiac tissue in contact with the balloon surface. It is used together with a compatible Medtronic 12 Fr inner diameter sheath (the sheath), the CryoConsole, and related components. For device compatibility questions, contact Medtronic Technical Support.

The balloon reaches cryoablation temperatures when refrigerant is injected from the CryoConsole to the balloon segment. A thermocouple positioned inside the balloon provides temperature reading capability. The catheter is introduced into the vasculature by traditional, minimally invasive techniques.

There are two radiopaque markers on the catheter to confirm the position of the balloon using fluoroscopy. The proximal radiopaque marker is located approximately 10 mm (0.394 in) proximal to the balloon. The distal radiopaque marker is located at the end of the injection tube.

Note: The 12 Fr FlexCath sheath is compatible with the catheter. There is one radiopaque marker on the FlexCath sheath, located approximately 5 mm (0.197 in) proximal to the tip of the sheath (see Figure 1). When the catheter's proximal radiopaque marker and the FlexCath sheath's radiopaque marker are aligned, the balloon is approximately 5 mm (0.197 in) outside of the FlexCath sheath. There are two shaft markers on the proximal section of the catheter shaft to visually confirm the position of the balloon within the FlexCath sheath. When the distal marker on the shaft of the catheter is aligned with the handle of the FlexCath sheath, the balloon segment is located inside the FlexCath sheath. When the proximal marker on the shaft of the catheter is aligned with the handle of the FlexCath sheath, the balloon segment is outside the FlexCath sheath (see Figure 1).

Figure 1. Arctic Front Advance cardiac cryoablation catheter

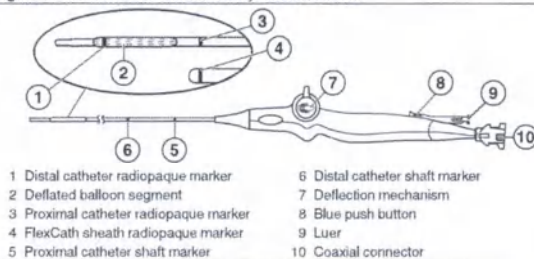The Arctic Front Advance cardiac cryoablation catheter is available in 2 models, as described in the following table:

Model	Inflated balloon diameter
2AF233	23 mm (0.91 in)
2AF283	28 mm (1.10 in)

For details about the CryoConsole and how to use it with the catheter to perform cryoablation procedures, see the *CryoConsole Operator's Manual*.

1.1 Contents of package

The catheter is supplied sterile. The package contains the following items:

- 1 Arctic Front Advance cardiac cryoablation catheter
- product documentation

2 Indications for use

The Arctic Front Advance cardiac cryoablation catheter is indicated for the treatment of atrial fibrillation (including paroxysmal, persistent, and long-standing persistent).

3 Contraindications

The Arctic Front Advance cardiac cryoablation catheter is contraindicated as follows:

- in the ventricle because of the danger of catheter entrapment in the chordae tendineae
- in patients with active systemic infections
- in patients with an intracardiac mural thrombus, myxoma, or interatrial baffle or patch
- in patients with cryoglobulinemia
- in patients with one or more pulmonary vein stents

4 Intended purpose

The intended purpose of the Arctic Front Advance cardiac cryoablation catheter is to ablate cardiac tissue in the left atrium.

5 Intended users

The intended user of the cryoablation system is a physician that performs cardiac ablation. Use of this device does not require specific qualifications or training beyond what is required per the applicable local, national, or federal regulations.

6 Intended patient population

The intended patient population of the Arctic Front Advance cardiac cryoablation catheter includes adults (patients 18 years or older) diagnosed with atrial fibrillation who are candidates for ablation. Ablation during pregnancy is not recommended. For the most up-to-date medical consensus for additional considerations, please refer to applicable treatment guidelines on the treatment of atrial fibrillation.

7 Clinical benefit

The expected clinical benefit to the patient includes restoring sinus rhythm, eliminating or reducing symptoms of atrial fibrillation (such as dizziness, palpitations, shortness of breath, rapid heartbeat or fatigue), and improving quality of life.

8 Summary of safety and clinical performance

The Summary of Safety and Clinical Performance (SSCP) can be found at <https://ec.europa.eu/tools/eudamed>. Search for the SSCP using the manufacturer and device name, and any of the following elements, as applicable: device model, reference number, catalog number, or the Basic Unique Device Identification (Basic UDI-DI) number — 0763000B00000256W.

9 Warnings and precautions

Anticoagulation therapy – Administer appropriate levels of peri-procedural anticoagulation therapy for patients undergoing left-sided and transseptal cardiac procedures. Administer anticoagulation therapy during and post-procedure according to the institution's standards. The Arctic Front Advance Cryoballoon was not studied for the safety of changes in anticoagulation therapy in patients with atrial fibrillation.

Balloon inflation or deflation – Do not inflate the balloon inside the sheath. Always verify with fluoroscopy or other appropriate visualization techniques that the balloon is fully outside the sheath before inflation to avoid catheter damage.

- Do not inflate the balloon while the catheter is positioned inside a pulmonary vein. Always inflate the balloon in the atrium and then position it at the pulmonary vein ostium. Inflating the balloon in the pulmonary vein may result in vascular injury.
- If the balloon cannot be inflated or deflated using the CryoConsole, have a Manual Retraction Kit on hand during the procedure. (See Chapter 13 for detailed instructions on the Manual Retraction Kit).

Biohazard disposal – Discard all used catheters and sterile components in accordance with hospital procedures.

Cardioversion or defibrillation during ablation procedure – Disconnect the catheter's electrical connection before cardioversion or defibrillation. Failure to do so may trigger system messages indicating a need for catheter exchange.

Catheter handling – Use extreme care when manipulating the catheter. Lack of careful attention may result in injury such as perforation or tamponade.

- Do not use excessive force to advance, withdraw, or rotate the catheter, especially if resistance is encountered. Excessive force may lead to catheter damage, including kinking of the guide wire lumen within the balloon segment.
- Do not use the catheter if it is kinked, damaged, or cannot be straightened.
- Straighten the shaft before inserting or withdrawing the catheter.
- Do not at any time preshape or bend the catheter shaft or balloon segment. Bending or kinking the catheter shaft may damage internal structures and increase the risk of catheter failure. Prebending of the distal curve may damage the catheter.

Note: Carefully grip the catheter shaft near the sleeve in order to help avoid kinking the catheter.

- Catheter advancement should be performed using fluoroscopy or other appropriate techniques.
- Do not position the cryoballoon catheter within the tubular portion of the pulmonary vein to minimize phrenic nerve injury and pulmonary vein stenosis.

Catheter integrity – Do not use the catheter if it is kinked or damaged. If the catheter becomes kinked or damaged while in the patient, remove it and use a new catheter. Before injecting, the physician should ensure that there is no kink in the catheter.

Circular mapping catheter compatibility – Use only Medtronic circular mapping catheters compatible with the inner lumen of the Arctic Front Advance Cryoballoon. Use of another mapping catheter may damage the catheter or compromise the procedure.

Contrast media – Use appropriate levels of contrast media in patients with comorbidities such as recent history of renal disease. Follow contrast labeling and institutional procedures regarding the appropriate medical strategies to minimize risk when using contrast media.

Correct guide wire or circular mapping catheter insertion and positioning – Do not advance the balloon beyond the guide wire or circular mapping catheter to reduce the risk of tissue damage.

- Ensure that the guide wire or circular mapping catheter is inserted into the catheter and through the balloon portion for adequate support during vascular access insertion. Failure to do so may result in catheter damage.

Cryoablation near prosthetic heart valves – Do not pass the catheter through a prosthetic heart valve (mechanical or tissue). The catheter may become trapped in the valve, damaging the valve and causing valvular insufficiency or premature failure of the prosthetic valve.

Cryo-adhesion – Do not pull on the balloon catheter, circular mapping catheter, sheath, umbilical cables, or CryoConsole while the balloon catheter or circular mapping catheter are frozen to tissue. This may lead to tissue injury. Before moving these components, use appropriate techniques to ensure that the balloon catheter and circular mapping catheter are not adhered to tissue.

Damage to lung or tracheobronchial tree – Damage to the lung or tracheobronchial tree has been observed in some subjects who have undergone left atrial ablation with the Arctic Front family. The physician should consider appropriate medical strategies to minimize the risk of damage to the lung or tracheobronchial tree.

Do not resterilize – Do not resterilize this device for the purpose of reuse. Resterilization may compromise the structural integrity of the device or create a risk of contamination of the device that could result in patient injury, illness, or death.

Embolism risk – Introducing any catheter into the circulatory system entails the risk of air or gas embolism, which may occlude vessels and lead to tissue infarction with serious consequences. Always advance and withdraw components slowly to minimize the vacuum created and therefore minimize the risk of air embolism.

Environmental limits – Perform cryoablation procedures only within the environmental parameters. Operating outside these parameters may prevent the start or completion of a cryoablation procedure. Refer to *Chapter 12, Specifications, page 8* for environmental parameters.

Esophageal injury – Esophageal ulcerations have been observed in some subjects who have undergone left atrial ablation with the Arctic Front family. As with other forms of left atrial ablation, the physician should consider appropriate medical strategies to minimize the risk of esophageal injury.

Fluid incursion – Do not expose the catheter handle or coaxial and electrical connectors to fluids or solvents. If these components get wet, the system may not function properly. This may lead to patient injury.

Fluoroscopy required for catheter placement – The use of fluoroscopy during catheter ablation procedures presents the potential for significant x-ray exposure to both patients and laboratory staff. Extensive exposure may result in acute radiation injury and increased risk for somatic and genetic effects. Only perform catheter ablation after giving adequate attention to the potential radiation exposure associated with the procedure, and taking steps to minimize this exposure. Give careful consideration before using the device in pregnant women.

For single use only – This device is intended only to be used once for a single patient. Do not reprocess or resterilize this device for the purpose of reuse. Reuse, reprocessing, or resterilization may compromise the structural integrity of the device or create a risk of contamination of the device that could result in patient injury, illness, or death.

Frequent flushing of the guide wire lumen – Flush the guide wire lumen before initial insertion and then frequently throughout the procedure to prevent coagulation of blood in the lumen. Flush the guide wire lumen with saline after each contrast injection.

Guide wire compatibility – Use only 0.81 mm (0.032 in) or 0.89 mm (0.035 in) guide wires with the catheter. Using another guide wire may damage the catheter or compromise the procedure.

Improper connection – Do not connect the cryoablation catheter to a radiofrequency (RF) generator or use it to deliver RF energy. This may cause catheter malfunction or patient injury.

Induced arrhythmias – Catheter procedures may mechanically induce arrhythmias.

Leakage current from connected devices – Use only isolated equipment (IEC 60601-1 Type CF equipment, or equivalent) with the CryoConsole and catheters or patient injury or death may occur.

Other catheters, devices, or wires – Avoid catheter entanglement with other catheters, devices, or wires. Such entanglement may necessitate surgical intervention.

Phrenic nerve impairment – To reduce the potential for phrenic nerve impairment, perform the following steps:

- Position the balloon as antral as possible and not in the tubular portion of the pulmonary vein. To ensure proper catheter position, observe the balloon shape and the balloon position within the anatomy.
- Monitor the status of the phrenic nerve continuously during right-sided pulmonary vein applications using an appropriate monitoring technique. One common monitoring technique is to continuously pace the phrenic nerve throughout each cryoablation application of the right pulmonary veins. While pacing, monitor diaphragm contraction by placing a hand on the abdomen to assess for loss of capture or changes in the strength of the diaphragmatic contraction.
- Stop ablation immediately if phrenic nerve impairment is observed.

Post-ablation period – Closely monitor patients undergoing cardiac ablation procedures during the post-ablation period for clinical adverse events.

Pressurized refrigerant – The catheter contains pressurized refrigerant during operation. Release of this gas into the circulatory system due to equipment failure or misuse may result in gas embolism.

Pulmonary vein narrowing or stenosis – Catheter ablation procedures inside or near pulmonary veins may induce pulmonary vein narrowing or stenosis. Do not ablate in the tubular

portion of the pulmonary vein. The occurrence of this complication may necessitate percutaneous angioplasty or surgical intervention.

Qualified users – This cryoablation system should be used only by physicians that perform cardiac ablation.

Required use environment – Cryoablation procedures should be performed only in a fully equipped facility.

RF ablation – Before powering up an RF generator or applying RF energy, disconnect the cryoablation catheter from the CryoConsole to avoid an error message and unnecessary catheter replacement.

Safety in magnetic resonance imaging (MRI) not evaluated – The Arctic Front Advance catheter has not been evaluated for safety and compatibility in the magnetic resonance (MR) environment. The safety of the Arctic Front Advance catheter in the MR environment is unknown (such as any heating, migration, or image artifact). Scanning a patient during the use of this device may result in patient injury.

Septal damage – Always deflate the balloon and withdraw the balloon into the transseptal sheath before removing the balloon from the left atrium. Crossing the septum while the balloon is unsheathed, inflated, or inflating in the septal puncture site may cause serious septal damage.

Serious incident – If a serious incident related to the device occurs, immediately report the incident to Medtronic and the applicable competent authority or regulatory body.

Steerable sheath compatibility – Use only a compatible Medtronic 12 Fr inner diameter sheath with the Arctic Front Advance Cryoballoon. Using another sheath may damage the catheter or balloon segment.

Sterile package inspection – Visually inspect all sterile-barrier packaging before use. If the device is damaged or the integrity of the sterilization barrier has been compromised, do not use the product. Contact your Medtronic representative for return information.

Storage and transport – There is no special handling, including no special storage or transport conditions, required for this device. Standard storage conditions are sufficient to safeguard the device. Store the device in the original packaging at room temperature in a dry place.

System compatibility – The Medtronic cryoablation catheters, refrigerant tanks, and other Medtronic CryoConsole components should only be used with the Medtronic CryoConsole.

Technical manual information – If you find information in this manual that is incorrect or illegible, contact your Medtronic representative or your local competent medical authority.

10 Adverse events

Potential adverse events associated with cardiac catheter cryoablation procedures include, but are not limited to, the following conditions:

- Access site complications (such as bruising or ecchymosis)
- Anemia
- Anxiety
- Arrhythmia (such as atrial flutter, bradycardia, heart block, or tachycardia)
- Back pain
- Bleeding from puncture sites
- Bronchial constriction
- Bronchial fistula
- Bronchitis
- Bruising
- Cardiac ischemia
- Cardiac tamponade
- Cardiopulmonary arrest
- Cerebral vascular accident
- Chest discomfort, pain, or pressure
- Cold feeling
- Coronary artery spasm
- Cough
- Death
- Diarrhea
- Dizziness
- Embolism
- Entrapment
- Esophageal damage (including atrio-esophageal fistula)
- Fatigue
- Fever
- Headache
- Hemoptysis
- Hypotension or hypertension
- Infection (for example, sepsis or urinary)
- Lightheadedness
- Nausea or vomiting
- Perforation
- Pericardial effusion
- Pericarditis
- Phrenic nerve injury
- Pleural effusion
- Pneumonia
- Pneumothorax
- Pseudoaneurysm
- Pulmonary edema
- Pulmonary hemorrhage
- Pulmonary vein dissection
- Pulmonary vein stenosis
- Renal insufficiency or renal failure
- Shivering
- Shortness of breath
- Sore throat
- Tissue infarction (such as myocardial infarction or renal infarction)
- Transient ischemic attack
- Vagal nerve injury (such as gastroparesis)
- Vasovagal reaction
- Visual changes (such as blurred vision)

11 Instructions for use

11.1 Connecting the Arctic Front Advance Cardiac Cryoablation Catheter

During connection of the catheter, keep the connections dry at all times. Fluid incursion may lead to system malfunction. To connect the catheter, follow these steps. (For more detailed instructions, see the *CryoConsole Operator's Manual*.)

1. Connect the auto connection box to the CryoConsole.
Note: The ECG cable is not required for an Arctic Front Advance Cryoballoon procedure, and should not be connected to the auto connection box.
2. Connect the catheter to a sterile electrical umbilical cable and then pass the other end of the cable out of the sterile field and connect it to the auto connection box.
3. Connect the catheter to a sterile coaxial umbilical cable and then pass the other end of the cable out of the sterile field and connect it to the CryoConsole.
4. Enable the vacuum on the CryoConsole.

11.2 Cryoablation

To use the catheter for a cryoablation procedure, follow these steps. (For more detailed instructions, see the *CryoConsole Operator's Manual*.)

Notes:

- Before introducing the catheter into the patient, test the deflection mechanism on the handle to ensure it is operational.
 - Always use the deflection mechanism on the handle to straighten the shaft before insertion or withdrawal of the catheter.
 - The Achieve family of mapping catheters is compatible for use with the Medtronic Arctic Front Advance Cryoablation Catheter and may be used to support and position the catheter.
1. Using aseptic technique, create a vascular access with an appropriate introducer. Obtain left atrial transseptal access using a transseptal sheath, its dilator, and needle.
 - Place the standard diagnostic pacing catheters.
 - Visualize the left atrial anatomy to help select a balloon size. Select the balloon size.
 - Selection of the balloon size should be based on pulmonary vein diameter and shape, the surrounding anatomy, and the desired position of the balloon outside the tubular portion of the pulmonary vein.
 2. Remove the transseptal sheath and dilator, leaving the guide wire positioned preferably in the left superior pulmonary vein.

3. Advance the sheath and dilator over the wire into the left atrium.
4. Slowly remove the guide wire and dilator from the sheath.
5. Aspirate and flush the sheath.
6. Obtain a catheter. Connect the Tuohy-Borst adapter and manifold to the push button Luer.
7. Flush the guide wire lumen, ensuring there are no air bubbles in the guide wire lumen or connections.
8. Load the circular mapping catheter or guide wire into the guide wire lumen.
9. Pull the circular mapping catheter or guide wire back into the balloon catheter until it is just inside the balloon catheter tip.
10. Flush the guide wire lumen again, ensuring there are no air bubbles in the guide wire lumen or connections.
11. Retract the sleeve from the balloon onto the catheter shaft.
12. Submerge and rinse the balloon to remove trapped air bubbles.
13. Advance the sleeve back over the balloon while the balloon is still submerged.
14. Place the sleeve adjacent to the hemostasis valve of the sheath.

Note: To avoid damaging the valve or introducing air, do not push the sleeve through the valve opening.

15. Insert the catheter into the sheath.

Notes:

- Carefully grip the catheter shaft near the sleeve to help avoid kinking the catheter.
- Follow institutional procedures to aspirate and flush the sheath.

16. Advance the balloon until the distal tip of the catheter aligns with the distal tip of the sheath using fluoroscopic guidance or other appropriate visualization techniques.
17. Advance the circular mapping catheter or guide wire to the target pulmonary vein using fluoroscopic guidance or other appropriate visualization techniques.
18. Advance the balloon over the circular mapping catheter or guide wire into the left atrium using fluoroscopic guidance or other appropriate visualization techniques.
19. Inflate the balloon in the left atrium by pressing the Start button on the CryoConsole control panel for 2 seconds.
20. Position the catheter at the ostium of the target pulmonary vein and not inside the tubular portion of the pulmonary vein.
21. Verify the balloon position by injecting a mixture of 50/50 contrast and saline into the catheter guide wire lumen port or by using other appropriate techniques. Be sure to flush the guide wire lumen with saline after each contrast injection.

Notes:

- To improve balloon position and support, reposition the circular mapping catheter. If a stable balloon position cannot be obtained, exchange the circular mapping catheter for a guide wire.
- Before exchanging the circular mapping catheter for a guide wire, retract the catheter into the sheath.
- When using an auto injector for contrast delivery, ensure that the pressure limit does not exceed 500 psig.
- Follow contrast labeling and institutional procedures regarding the appropriate medical strategies to minimize the risk to the patient associated with using contrast.

22. Perform the cryoablation.

Notes:

- During the initial ablation phase, monitor the balloon's position using appropriate visualization techniques. Moisture in the system may cause the inflated balloon to deflate and then re-inflate during the initial transition phases.
- Set the ablation duration on the CryoConsole screen.
- Physicians may modify the preset ablation duration based on clinical judgment. It is not recommended that the ablation duration exceeds 240 seconds.
- The balloon's outer diameter varies from inflation to cryoablation and may cause the balloon to shift.

23. Wait for the cryoablation phase to complete (at the end of the preset duration). The balloon remains inflated and the thawing phase begins.

Note: At any time, the ablation can be stopped by pressing the Stop Current Action button on the CryoConsole control panel or by pressing the Stop CryoAblation button on the screen.

24. During the thawing phase, observe the temperature indicator on the screen. The balloon deflates automatically when the temperature reaches 20°C.

25. Before moving the balloon, use appropriate techniques to ensure that the catheter is not adhered to tissue.

26. Determine effective ablation of the cardiac tissue by assessing electrical isolation of the pulmonary vein from the left atrium (entrance and exit block) after the cryoablation is complete.

Note: As needed, perform additional treatments by positioning the balloon differently in the same pulmonary vein.

27. Position the catheter at the ostium of the next target pulmonary vein using the sheath, circular mapping catheter, or guide wire. Return to Step 17 and continue ablation.

28. To retract the balloon into the sheath, perform the following steps:

Note: Ensure that the distal tip of the catheter is free to move to its maximum length.

- a. Use the deflection mechanism on the sheath handle to straighten the sheath. Use the deflection mechanism on the catheter handle to straighten the catheter.

- b. Inflate the balloon.

- c. Perform the following two steps simultaneously:

Advance the blue push button on the catheter handle, as shown in Figure 2. This causes the balloon to extend to maximum length and wrap tightly.

Deflate the balloon by pressing the Stop Current Action button on the CryoConsole control panel, or by selecting the Deflates Balloon option on the screen.

- d. Retract the catheter into the sheath.

Figure 2. Catheter handle

29. Remove the catheter from the patient.

30. After the procedure, follow the instructions to "Shut down the system" in the *CryoConsole Operator's Manual*. Follow the prompts to close the refrigerant tank and replace the cap on the coaxial port. Confirm these actions when prompted by the CryoConsole.

12 Specifications

Significant physical and performance characteristics are included in the table below.

Description	Specification
Catheter shaft outer diameter	3.5 mm (10.5 Fr; 0.140 in)
Tip outer diameter	3.0 mm (9 Fr; 0.118 in)
Tip length	13.5 mm (0.530 in)

Guide wire lumen inner diameter	1.27 mm (0.050 in) nominal
Inflated balloon diameter:	
2AF233	23 mm (0.91 in) nominal
2AF283	28 mm (1.10 in) nominal
Effective catheter length (with balloon inflated)	95.0 ± 2.0 cm (37.40 ± 0.80 in)
Number of thermocouples	1
Compatibility:	
Console	Medtronic CryoConsole*
Electrical connection cable	Medtronic electrical umbilical cable (2035U, 2035UC)
Coaxial connection cable	Medtronic coaxial umbilical cable (203CX, 203CXC)
Introducer sheath	Medtronic 4.0 mm (12 Fr; 0.158 in) inner diameter
Guide wire	0.81 mm (0.032 in) and 0.89 mm (0.035 in)
Mapping catheter	Medtronic circular mapping catheters with a maximum outer diameter of 1.09 mm (0.043 in)
Operating environmental parameters	15°C to 30°C (59°F to 86°F) at altitudes less than 2400 m (8000 ft) above sea level at humidity between 30% and 75%
Sterilization	Single use, sterile device, ethylene oxide

* The specific model number for the compatible CryoConsole may vary based on geography. Contact your local Medtronic representative for additional information.

13 Accessories

13.1 Manual Retraction Kit


13.1.1 Description

The Model 20MRK manual retraction kit (MRK) is comprised of a male coaxial connector, a stopcock, and a vacuum-lock syringe (Figure 3). It is a sterile device intended for single use. The MRK connects to the Arctic Front family of cardiac cryoablation catheters. It provides an alternative means of inflating and deflating the catheters, for balloon wrapping and retraction, in unlikely situations when standard catheter retraction methods cannot be used.

The MRK is designed for use with the Arctic Front family of cardiac cryoablation catheters.

Note: The MRK is only used with the Arctic Front family of cardiac cryoablation catheters. Please refer to the catheter information related to the indications for use, contraindications, intended users, adverse events, warnings and precautions, intended patient population, and clinical benefit.

Figure 3. Manual Retraction Kit

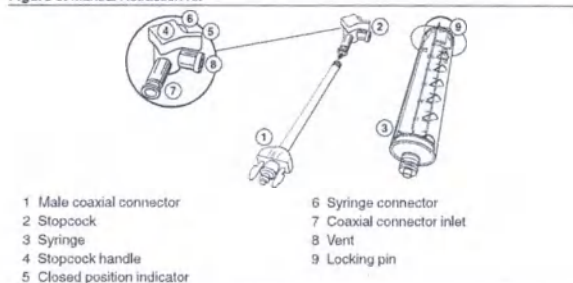The manual retraction kit has the following parts:

- **Male coaxial connector** – connects the manual retraction kit to the female coaxial connector of the catheter.
- **Stopcock** – provides the user an ability to control the direction of flow into the catheter or to the atmosphere.
- **Syringe** – provides the user an ability to manually inflate or deflate the balloon by moving the plunger. Retracting the plunger creates a vacuum in the balloon and the catheter return path, resulting in balloon deflation. Balloon inflation occurs when saline is injected into the balloon and catheter return path.

13.1.2 Intended purpose

The intended purpose of the MRK is to aid manual retraction of the balloon catheter in unlikely situations when standard catheter retraction methods cannot be used.

13.1.3 Instructions for use

Notes:

- Two people are required for this procedure.
- Do not disconnect the coaxial and electrical umbilical cables until indicated in these steps.


13.1.3.1 Prepare the equipment

1. Unpack and assemble the manual retraction kit.
2. Verify that the stopcock and syringe are assembled as shown.


13.1.3.2 Release the vacuum and attach the coaxial connector

1. On the Therapy tab of the CryoConsole, touch **Release Vacuum** on the screen.
2. Disconnect the electrical umbilical cable and the coaxial umbilical cable from the catheter, as shown:

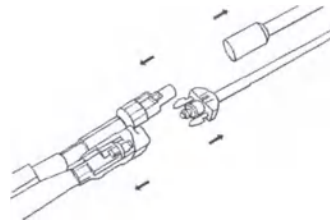3. Push the syringe plunger in.
4. Connect the manual retraction kit coaxial connector to the catheter coaxial connector:

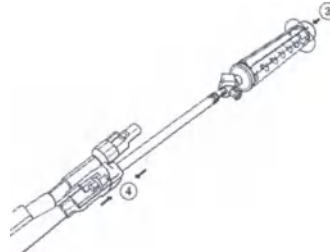13.1.3.3 Apply vacuum to separate the catheter

1. Rotate the stopcock handle to close flow to the vent:

2. Pull the plunger out of the syringe to approximately 60 cm³.
3. Rotate the plunger to lock the syringe:

4. Rotate the stopcock handle to close flow to the catheter inlet:


5. Disconnect the syringe from the stopcock.
6. Rotate the plunger to unlock it.
7. Push in the syringe plunger to empty the syringe.
8. Reconnect the syringe to the stopcock.
9. Apply a second syringe stroke of vacuum to the catheter by repeating Step 1 through Step 8.

13.1.3.4 Prepare and perform saline inflation of the balloon

1. Ensure the flow to the catheter inlet is closed:

2. Disconnect the syringe from the stopcock:

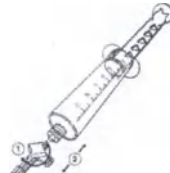3. Fill the syringe with 30-50 cm³ of saline.
Note: Minimum quantity is 20 cm³.


4. Connect the syringe to the stopcock:

5. Rotate the stopcock handle to close flow to the vent:


6. Inject the saline into the catheter (minimum 20 cm³):

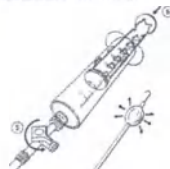13.1.3.5 Tension and deflate the balloon

1. Retract the saline, using two people, as follows:
 - a. Advance the blue push button on the catheter handle.
 - b. At the same time, pull on the plunger to retract the saline.

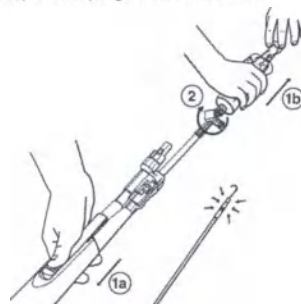2. Rotate the stopcock handle to close flow to the catheter inlet:

3. Follow steps 28 and 29 in Section 11.2, Cryoablation, page 7 to remove the catheter from the patient.

Caution: Discard all used catheters and sterile components in accordance with hospital procedures.

13.1.4 Cleaning and Storage

Before use, store the manual retraction kit in accordance with temperature guidelines described in "Storage and transport", Chapter 9, Warnings and precautions, page 5.

The manual retraction kit is a sterile component intended for single use. Do not attempt to clean it.

13.1.5 Disposal

Do not dispose of this product in the unsorted municipal waste stream. Dispose of this product according to local regulations. See <http://recycling.medtronic.com> for instructions on proper disposal of this product.

Caution: Discard all used catheters and sterile components in accordance with hospital procedures.

13.1.6 Specifications

Table 1. Manual retraction kit specifications

Description	Specification
Components	Medtronic male coaxial connector 3-way, high pressure stopcock 60 cm ³ , vacuum lock syringe
Compatibility	Medtronic Arctic Front family of cardiac cryoablation catheters
Operating environmental parameters	15°C to 30°C (59°F to 86°F) at altitudes less than 2400 m (8000 ft) above sea level at humidity between 30% and 75%
Sterilization	Single use, sterile device, ethylene oxide

14 Medtronic limited warranty

For complete warranty information, see the accompanying limited warranty document.

15 Service

Medtronic employs highly trained representatives and engineers located throughout the world to serve you and, upon request, to provide training to qualified hospital personnel in the use of Medtronic products. Medtronic also maintains a professional staff to provide technical consultation to product users. For more information, contact your local Medtronic representative, or call or write Medtronic at the appropriate telephone number or address listed on the back cover.

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© 2019 Medtronic
Medtronic, Medtronic with rising man logo, and Medtronic logo are trademarks of Medtronic. Third-party trademarks ("TM") belong to their respective owners. The following list includes trademarks or registered trademarks of a Medtronic entity in the United States and/or in other countries.
(Medtronic B.V.)