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User Manual

**Device for extracorporeal shock wave therapy (ESWT),
model SALUS - RSWT**

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The Instructions For Use provides a guide to the efficient use of the SALUS-RSWT. For clinical definitions, pathological effects etc., of its functions, refer to the related medical publications.

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If you encounter any problems during operation of the equipment, stop its use immediately, and contact the customer service center of REMED (Refer to chapter 14).

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Contents

1. Information on safety	7
1.1 Definitions of symbols	7
1.2 Requirements for safety	9
1.2.1 Precautions for electrical safety	9
1.2.2 Classifications	9
1.3 Precautions during operation	10
1.4 Indications for use, Contraindications for use, Possible complications and Possible side effects*	11
1.5 Abnormal reaction	11
1.6 General warnings	11
2. Extracorporeal Shockwave (Learn RSWT)	12
2.1 Mechanism	12
2.2 Objective of application	13
2.3 Intended use for RSWT	13
3. Installation	14
3.1 Components	14
3.2 Electrical installation condition	15
3.3 Environmental condition	15
3.3.1 Application environment	15
3.3.2 Transport and Storage environment	15
3.4 Installation method	15
3.4.1 Precautions during installation	15
3.4.2 Connection of power	16
3.4.3 Moving and fixing of equipment	16
3.4.4 Gun	17
4. Product description	18
4.1 External view	18
4.1.1 Front view	18
4.1.2 Rear view	18
4.2 Gun	19
4.3 Controls	19
5. Preparations for operation	20
5.1 Moving and fixing equipment	20
5.2 Power connection	20
6. Operation instructions	21
6.1 Applications	21
6.2 Operation mode	22

6.3 Environment setting	26
7. Messages	27
7.1 Over Temperature	27
7.2 Gun Open	27
7.3 Components exchange messages	28
8. Treatment protocols	29
9. Technical specifications	30
9.1 Size and weight	30
9.2 Power	30
9.3 Number of shocks	30
9.4 Modes	30
9.5 Intensity	30
9.6 Frequency	30
9.7 Applicator	30
10. EMC emissions and immunity	31
11. FAQ	37
12. Maintenance	38
12.1 Routine inspection of equipment	38
12.2 Safety inspection	38
12.3 Cleaning procedure	39
12.3.1 General cleaning	39
12.3.2 Cleaning applicators	39
12.3.3 Projectile & projectile guide cleaning	41
12.3.4 Connection part of the gun cleaning	41
12.4 Replacement of components	42
12.5 Removing air compressor condensate	43
12.6 Troubleshooting	44
12.7 Regular inspection for performance maintenance	46
12.8 Warranty	46
14. Contact information	48

List of figures

Figure 1. Mechanism of shockwave therapy	12
Figure 2. Main components	14
Figure 3. Power cable connection.....	16
Figure 4. Transporting wheel in 'Loose' (Left) and in 'Locked' (right)	17
Figure 5. Front view.....	18
Figure 6. Rear view	18
Figure 7. Gun.....	19
Figure 8. Controls.....	19
Figure 9. Application of the SALUS-RSWT	21
Figure 10. REMED logo	22
Figure 11. SALUS-RSWT initial screen.....	22
Figure 12. Main screen.....	22
Figure 13. Over temperature message.....	27
Figure 14. Gun open message	27
Figure 15. Tube & Projectile exchange message	28
Figure 16. Handpiece (Gun) exchange message	28
Figure 17. Applicator cleaning.....	39
Figure 18. 6 mm Applicator cleaning.....	40
Figure 19. 10 mm Applicator cleaning.....	40
Figure 20. 15 mm Applicator cleaning.....	40
Figure 21. 20 mm Applicator cleaning.....	40
Figure 22. 35 mm Applicator cleaning.....	41
Figure 23. Projectile & projectile guide cleaning	41
Figure 24. Gun connection part cleaning	41
Figure 25. Projectile, projectile guide, bumper	42
Figure 26. Projectile, projectile guide, bumper disassembling	42
Figure 27. Air compressor condensate remove.....	43
Figure 28. Hose connection	44
Figure 29. Hose disconnection	44

List of tables


Table 1. Treatment protocols	29
Table 2. Self-troubleshooting.....	44


1. Information on safety

1.1 Definitions of symbols

- ☐ When there is specific information in the Instructions For Use that needs to be emphasized for safety, the following terms and symbols are indicated. All warnings and precautions shall always be observed.
- ☐ The manufacturer or agent of the product is not responsible for any personal/material damage caused by erroneous use, operation for purposes other than its intended objective and negligence of product maintenance.

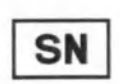	Warning	"Warning" symbol is used to indicate a life-threatening risk to the operator if the warning is disregarded.
	Caution	"Caution" symbol is used to indicate that injury or damage can be caused if the caution is disregarded.
	No pushing	To prohibit pushing against an object
	No sitting	To prohibit sitting on a surface
	BF Type	Applied part BF type
	Operation instructions	Operation instructions
	instructions for use.	Follow instructions for use.
	Serial number	Serial number
	Article number	Article number
	Manufacturer	Manufacturer
	Date of manufacture	Date of manufacture

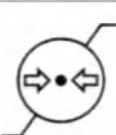	Do not stack	To indicate that the items shall not be vertically stacked, either because of the nature of the transport packaging or because of the nature of the items themselves.
	Keep away from rain	To indicate that the transport package shall be kept away from rain and in dry conditions.
	Fragile	To indicate that the contents of the transport package are fragile and the package shall be handled with care.
	This way up	To indicate correct upright position of the transport package.
	Recycling	To indicate the location of a recycling bin or container.
	"OFF" (Power)	To indicate disconnection from the mains, at least for mains switches or their positions, and all those cases where safety is involved.
	"ON" (Power)	To indicate connection to the mains, at least for mains switches or their positions, and all those cases where safety is involved.
	Alternating current	To indicate on the rating plate that the equipment is suitable for alternating current only to identify relevant terminals.
	WEEE	Indicates that when the end-user wishes to discard this product, it must be sent to separate collection facilities for recovery and recycling.
	Temperature limitation	To indicate the maximum and minimum temperature limits at which the item shall be stored, transported or used.
	Humidity Limitation	To indicate the acceptable upper and lower limits of relative humidity for transport and storage.
	Atmospheric pressure limitation	To indicate the acceptable upper and lower limits of atmospheric pressure for transport and storage.
	Dangerous Voltage	Warning, electricity

1.2 Requirements for safety

1.2.1 Precautions for electrical safety

- ☐ The rated power for this equipment is AC 220 V.
- ☐ Confirm that all the connecting sections (power line or peripheral equipment) are adequately connected to the equipment.
- ☐ Confirm that the equipment is completely grounded.
- ☐ Repair, expansion, and installation of equipment shall not be performed by anyone other than the specialized personnel authorized by the manufacturer. Arbitrary disassembling/assembling of equipment by the user is absolutely prohibited.
- ☐ Before connecting other equipment that is not specified in the Instructions For Use, be sure to inform this company or the agency with authority over product marketing.
- ☐ To avoid electrical noise during use, the equipment shall be installed at a significant distance from any generator, X-ray equipment, broadcasting device, mobile electric wire, and so forth.
- ☐ Needs attention due to potential electromagnetic or other interference between ME EQUIPMENT and other devices.
- ☐ An independent power circuit is essentially required, and sharing of power with other electronic devices is not recommended.
- ☐ Do not place the ME EQUIPMENT in a location where it is difficult to disconnect the power plug or other detachable plug.

1.2.2 Classifications

- ☐ Protection type and level for electric shock : Class I, Type BF
- ☐ Electromagnetic compatibility (EMC) test standard: Type 2, class A

1.3 Precautions during operation

- ☐ Installation and re-installation of equipment shall always be performed by specialized personnel authorized by the manufacturer.
- ☐ Installed equipment shall be subject to regular safety inspections from specialized personnel authorized by the manufacturer.
- ☐ Repair or installation of equipment shall only be performed by the specialized personnel authorized by the manufacturer, and as such, all arbitrary disassembly/assembly of equipment by the user is absolutely prohibited.
- ☐ For electrical safety, equipment shall always be connected to a safely grounded power supply for operation.
- ☐ The operator and manager of equipment shall thoroughly understand the Instructions For Use, and it must be kept in the close proximity of the equipment.
- ☐ To prevent safety accidents and ensure adequate maintenance, place the safety marks, guide phrases and regular inspection table provided with the equipment at an easily visible location near the equipment.
- ☐ As a fan is installed in the body of the equipment for air circulation, remove curtains or any other objects which can block the flow of air from near the equipment.
- ☐ From the surroundings of the location where equipment is installed, remove water, alcohol, flammable materials, and so forth.
- ☐ When equipment is operated, do not use any mobile phones, radio sets, mobile radio transmitters, wireless toys, etc. in the proximity of the equipment.
- ☐ During the operation of the equipment, the patient shall not take drinks, water, etc. which can influence the equipment.
- ☐ Do not separate gun and applicator during operation of the equipment, as it may cause injury due to recoil of air pressure.
- ☐ If you do not apply the Transmission Gel enough or do not come into contact with the skin, patients may feel pain because the air bubble is formed and the shock wave is generated. Therefore, it is necessary to apply gel sufficiently and close the gun to patients tightly to prevent side effects.
- ☐ Shockwave may be affected by overpressure when used around the lungs, ear, and gastrointestinal tract.
- ☐ Do not directly apply shockwaves to the heart, head or wounded area.
- ☐ This equipment should be used for patients who can accurately express and communicate pain during operation.

- ☐ Excessive pressure on the stimulation area with a gun may cause skin damage such as blisters and redness.
- ☐ Do not move the unit during use because the gun may fall or the unit may be damaged due to impact or movement of the unit.

1.4 Indications for use, Contraindications for use, Possible complications and Possible side effects.

- ☐ See the Appendix to the Instruction Manual for full details.

1.5 Abnormal reaction

- ☐ If any abnormality such as skin redness, pain, petechiae hematoma, hematoma, hemorrhage, edema, etc. is found, please stop use immediately and consult with your doctor.

1.6 General warnings

- ☐ Be sure to use the equipment in accordance with the Instructions For Use.
- ☐ In the Instructions For Use, the method for the correct use of the equipment is described. Read the Instructions For Use carefully.
- ☐ The equipment shall be operated under the supervision of a person who has completed specialized medical education.
- ☐ This company is not responsible for any disadvantages / damages related to the operation of equipment by a person without a medical education.
- ☐ This equipment shall not be altered, remodeled, or used for any purpose other than the intended objective.

2. Extracorporeal Shockwave (Learn RSWT)

2.1 Mechanism

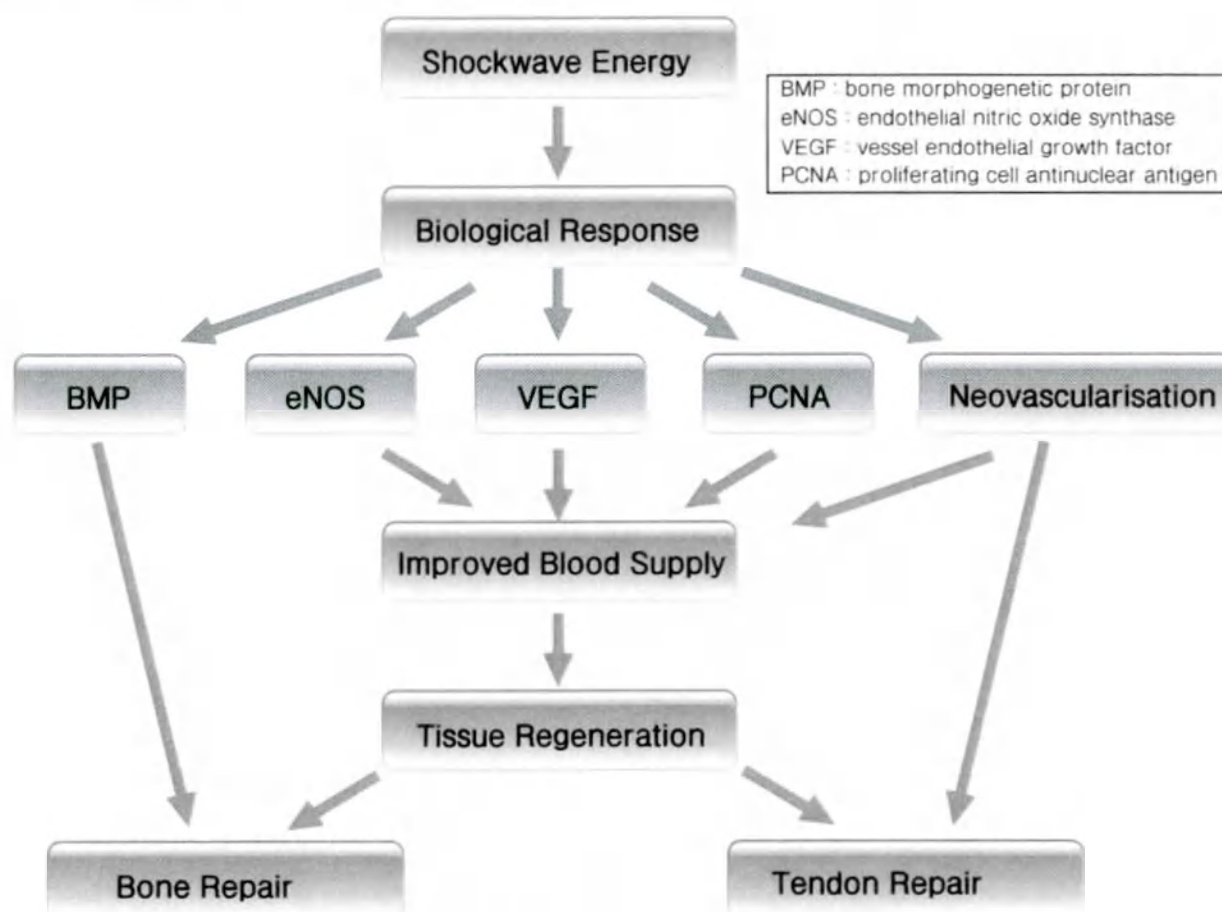
Figure 1. Mechanism of shockwave therapy

An extracorporeal shockwave is a waveform with strong pressure that occurs when the energy is suddenly released in a small space for a short period of time, which stimulates the damaged tissue and results in a therapeutic effect. Extracorporeal shockwaves stimulate the early expression of factors related to angiogenesis, eNOS (endothelial nitric oxide synthase, VEGF (vascular endothelial growth factor), and PCNA (proliferating cell antinuclear antigen). It promotes neovascularization and affects blood supply and cell proliferation. As a result, it activates surrounding tissues and bones to reduce pain and improve natural healing effects. It also helps to recover the bones.

The extracorporeal shockwave is also known to have a therapeutic effect on chronic inflammation. Chronic inflammation is caused by the fact that the inflammation reaction has not been completely interrupted, which can cause injury to healthy tissue, leading to chronic inflammation. Mast cells are the most important factor in the inflammation process. The extracorporeal shockwave activates these cells, stopping the inflammatory reaction, and helping with normal treatment and recovery processes. The extracorporeal shockwave also affects collagen which is the most important factor in the process of recovery of the upper musculoskeletal and ligament. It also promotes the formation of procollagen, making the

newly formed tendon more firm and firm. In addition, the extracorporeal shockwave is known to help treat pain by reducing the concentration of Substance P that is associated with pain signals.

2.2 Objective of application

This device is an extracorporeal shockwave device used for the pain relief. It is prohibited to use it for purposes other than the intended use.

2.3 Intended use for RSWT

The RSWT is a non-surgical alternative for the treatment of chronic proximal plantar fasciitis for patients 18 years of age or older with symptoms for 6 months or more and a history of unsuccessful conservative therapy.

3. Installation

3.1 Components

Main body with gun		Power cable		User guide		
				 <i>SALUS-RSWT</i> <i>Instructions For Use</i> <i>Extracorporeal Shock Wave Therapy</i>  <i>REMEDI Co., Ltd.</i>		
Gun		Applicator				
		6mm	10mm	15mm	20mm	35mm(option)
						
Cleaning kit						
Cleaning stick						
						
						
Tweezers						
			Hex wrench			

Figure 2. Main components

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3.4 Installation method

3.4.1 Precautions during installation

- ☐ Install the equipment on a flat surface.
- ☐ Confirm that the power is connected.
- ☐ Operate the SALUS-RSWT within the ambient temperature range of 10 – 30 °C, and a humidity between 30 – 85 %.
- ☐ Install the equipment in a location with appropriate ambient temperature and humidity, and do not install in a location where it will be exposed to dust or flammable materials.
- ☐ Be careful not to damage the equipment by excessive shock.

3.4.2 Connection of power

- ☐ Confirm that power socket and power cable are grounded before operating the machine.
- ☐ Plug the power cable into the power terminal of the SALUS-RSWT as shown in below, and connect the other end to a grounded power socket.
- ☐ Confirm that the power cable is properly connected to the power terminal and power socket. If the connection is incorrect, unexpected problems can occur.
- ☐ Do not connect several devices to one power outlet.
- ☐ In particular, do not use a power cable that generates connection noise.

Figure 3. Power cable connection

	Warning	To prevent the risk of electric shock, it must be connected to a protected grounded power source.
--	----------------	---

3.4.3 Moving and fixing of equipment

- ☐ Before moving the equipment, remove the connected power cable and the connection cables for peripheral equipment.
- ☐ When moving the equipment, place the transporting wheel at the lower section in the “loose” position.
- ☐ When the movement is completed, prevent the equipment from shaking by placing the transporting wheel in the “locked” position.
- ☐ If the equipment is moved while the wheel is in the “locked” position, it will damage the wheel, so please be careful.

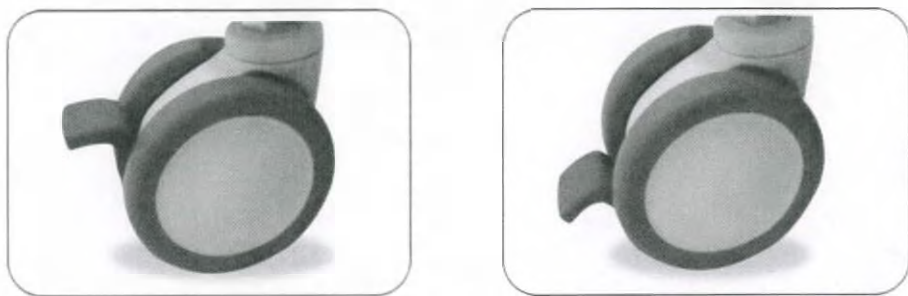
Figure 4. Transporting wheel in 'Loose' (Left) and in 'Locked' (right)

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	Warning Do not separate gun and applicator during the operation.
--	---

4. Product description

4.1 External view

4.1.1 Front view

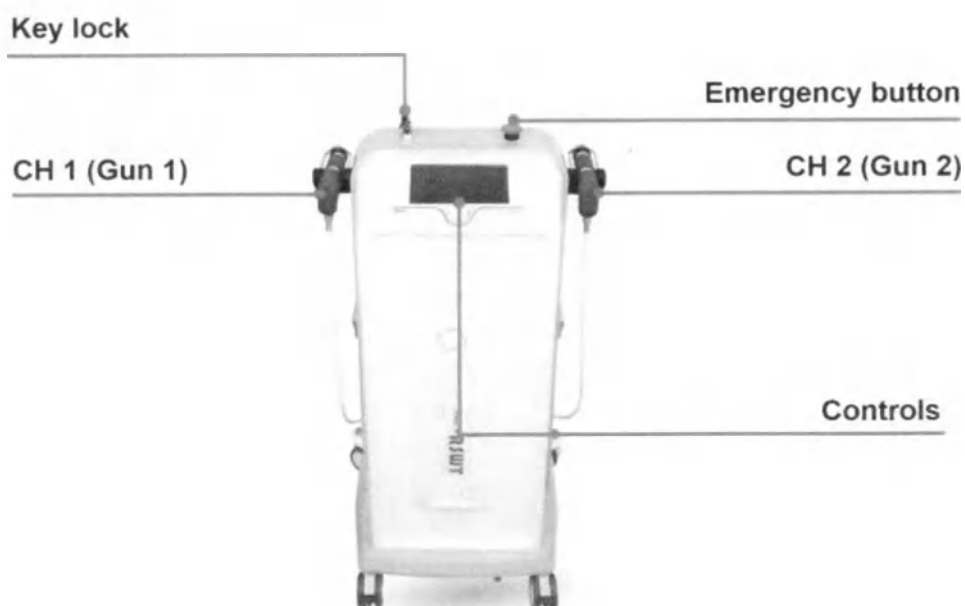
Figure 5. Front view

4.1.2 Rear view

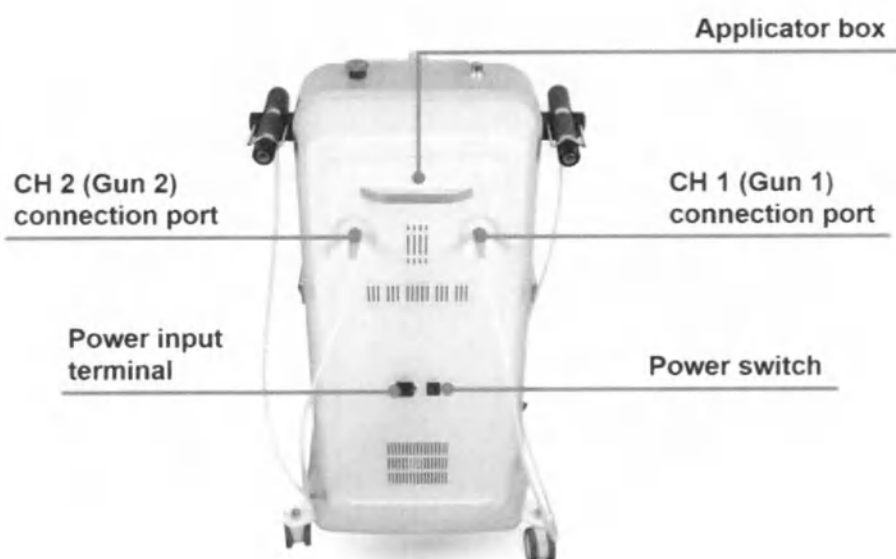
Figure 6. Rear view

4.2 Gun

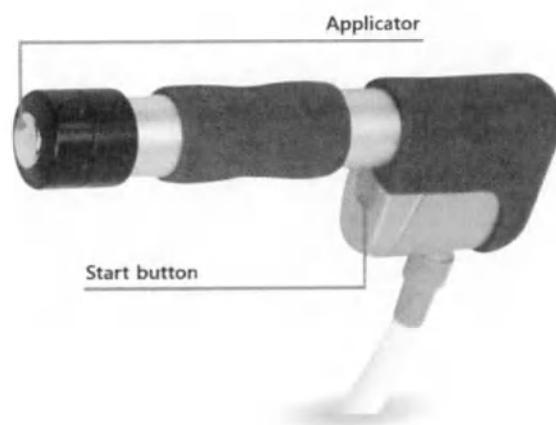
Figure 7. Gun

4.3 Controls

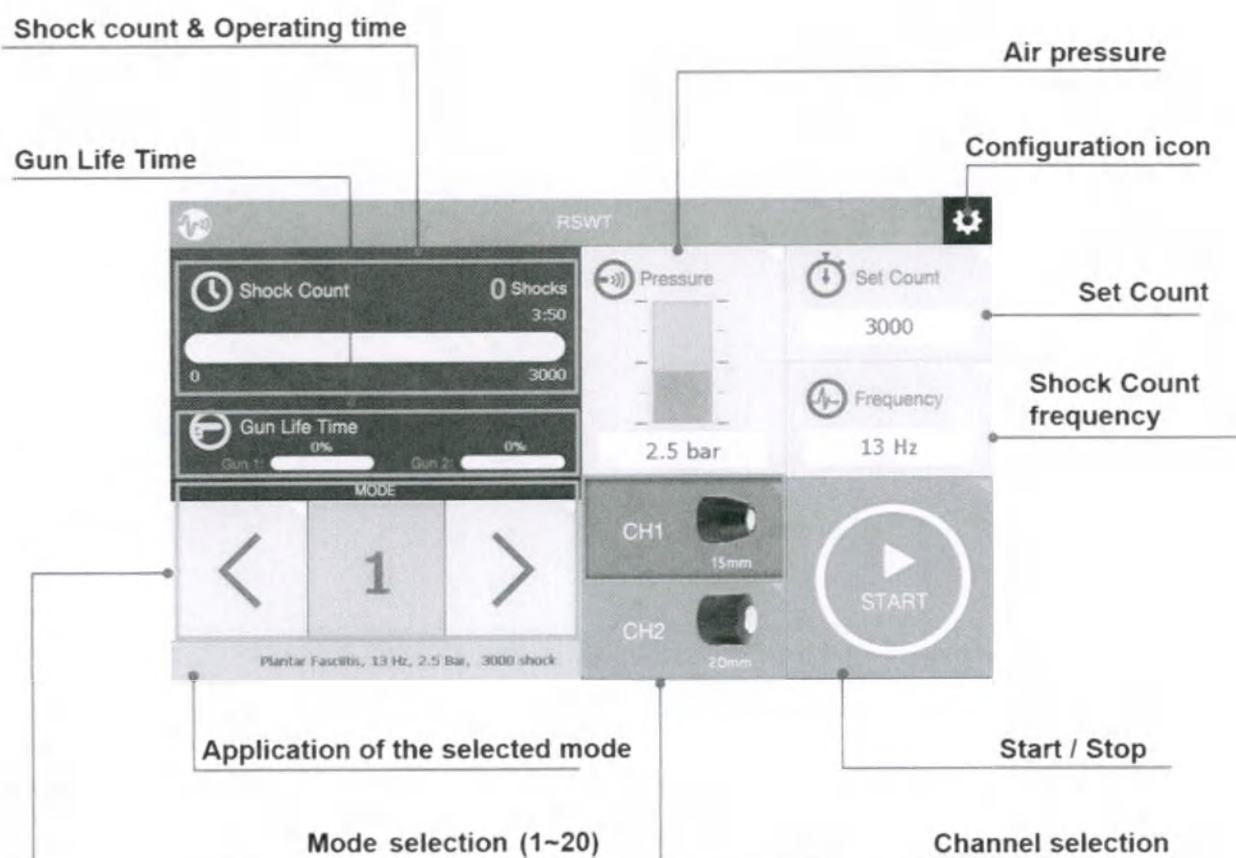
Figure 8. Controls

5. Preparations for operation

5.1 Moving and fixing equipment

- ☐ Before moving the equipment, remove the connected power cable and the connection cables for peripheral equipment.
- ☐ When moving the equipment, place the transporting wheel at the lower section in the "loose" position.
- ☐ When the movement is completed, prevent the equipment from shaking by placing the transporting wheel in the "locked" position.
- ☐ If the equipment is moved while the wheel is in the "locked" position, it will damage the wheel, so please be careful.

5.2 Power connection

- ☐ Confirm that power socket and power cable are grounded before operating the machine.
- ☐ Plug the power cable into the power terminal of the SALUS-RSWT, and connect the other end to a grounded power socket.
- ☐ Confirm that the power cable is properly connected to the power terminal and power socket. If the connection is incorrect, unexpected problems can occur.
- ☐ Do not connect several devices to one power outlet.
- ☐ In particular, do not use a power cable that generates connection noise.

6. Operation instructions

6.1 Applications

The SALUS-RSWT can be used at different sites of the body depending on the pain. Figure 9 shows the use of the most common applications.

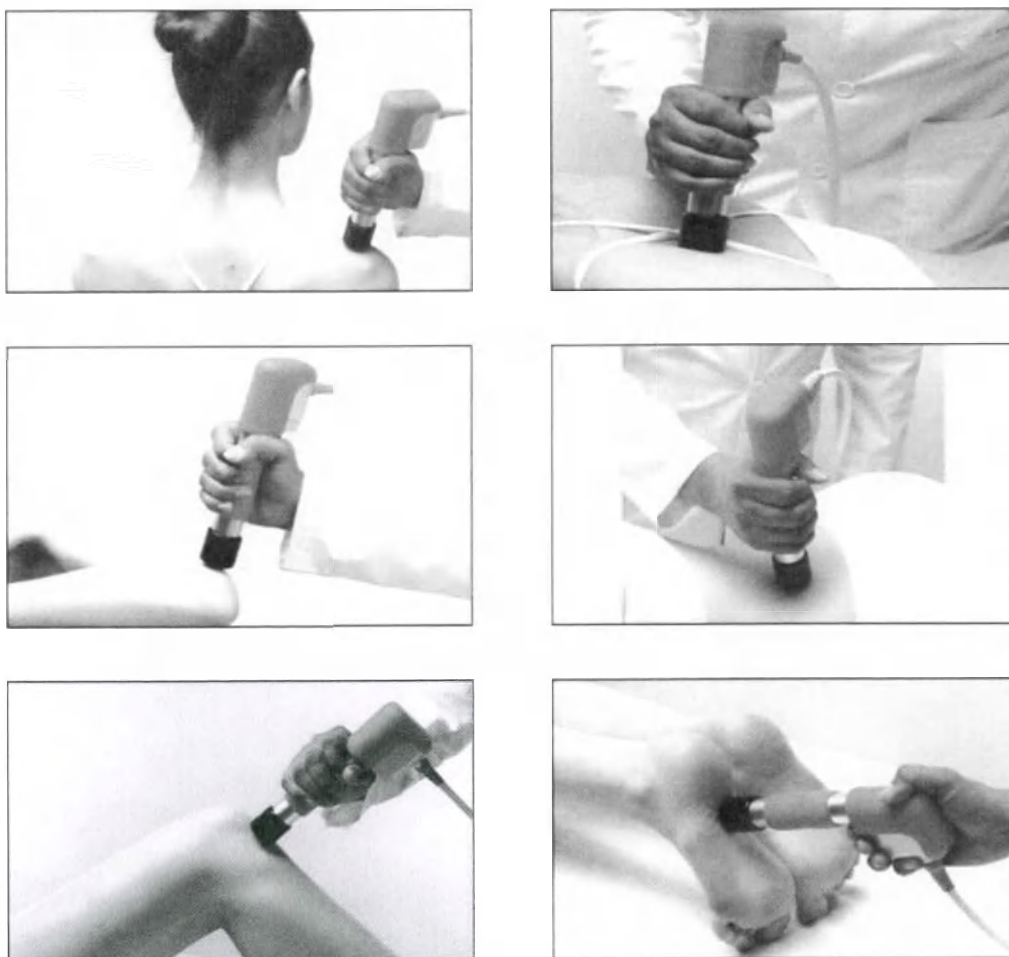
Figure 9. Application of the SALUS-RSWT

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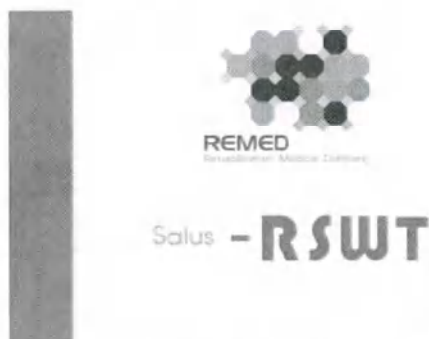
Figure 10. REMED logo

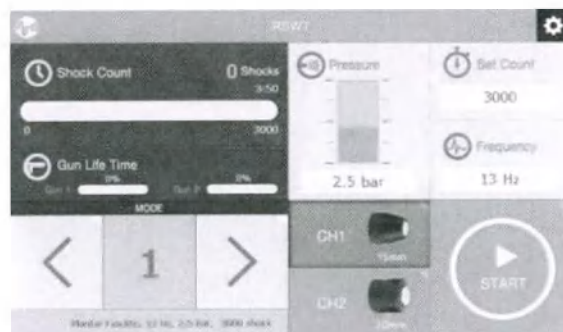
Figure 11. SALUS-RSWT initial screen

The compositions of the main screen are like the following.


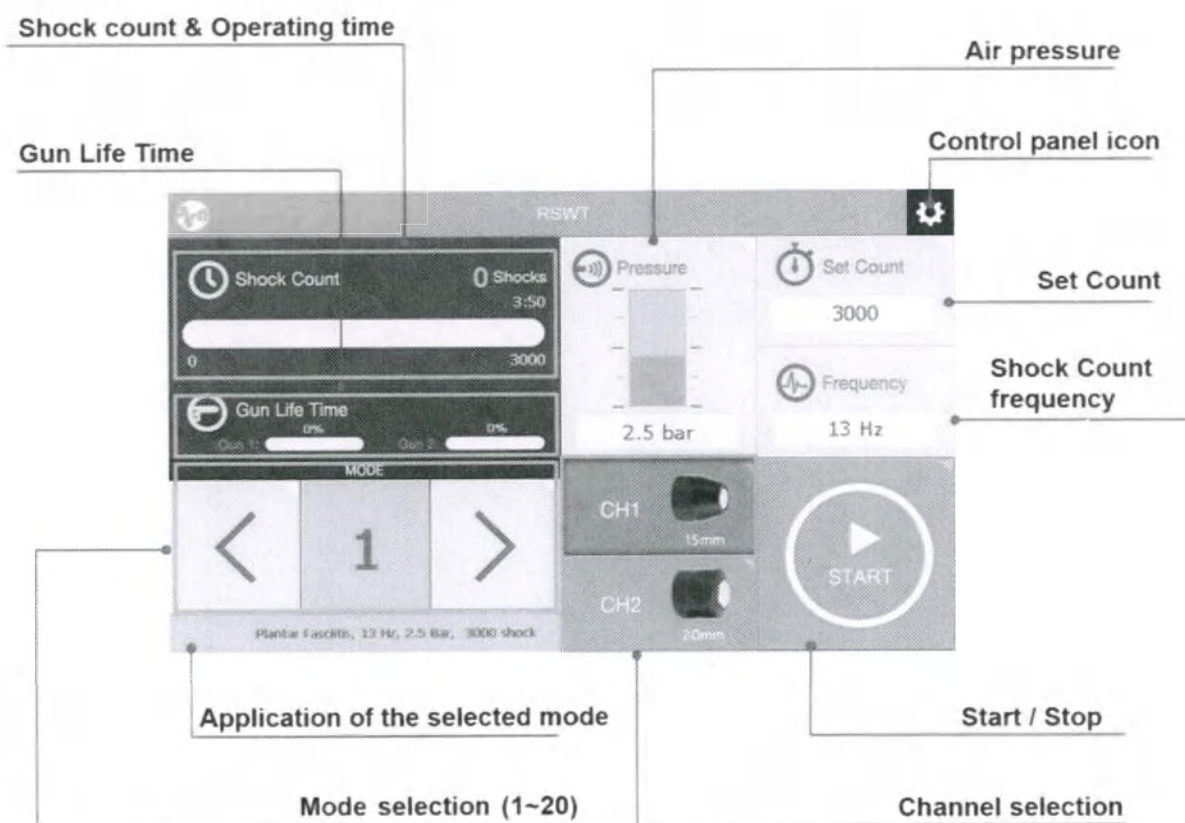
Figure 12. Main screen

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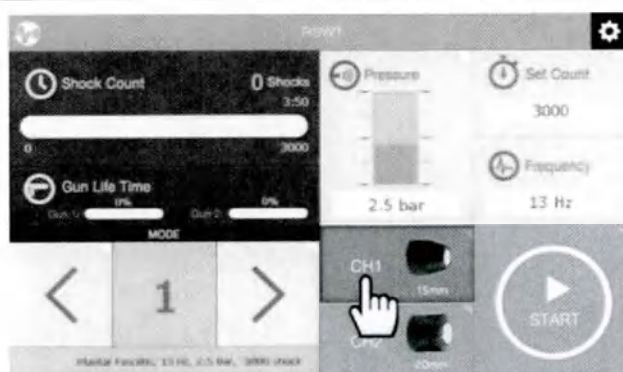
(a) Select the channel to use (The button of the activated channel becomes darker). Two channels can't be used simultaneously.

(b) Select the mode by using a left/right arrow or turning encoder. Application of the selected mode, parameter values and type of the applicator are shown.

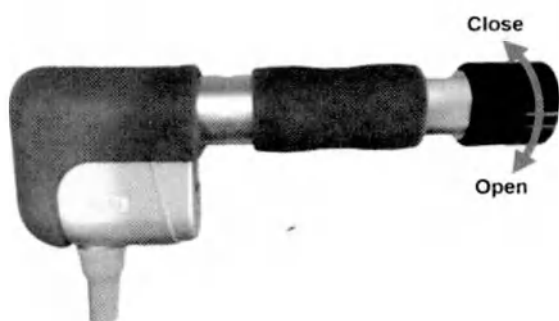
(c) Disconnect the applicator connected with the Gun by turning clockwise the applicator. Install the applicator selected in (b) by turning anti-clockwise the new applicator.

(d) To temporarily change the parameter value, press the window of each variable to select it and turn the encoder to change the value. However, when you turn the unit off and back on, or change to another channel, and then return to the original channel, the parameter value that was temporarily changed is not saved, so the value before the change appears. (E) to (j) are methods to save parameters after changing parameters


(e) When the user want to change the parameter values, press mode number for 2 ~ 3 seconds.

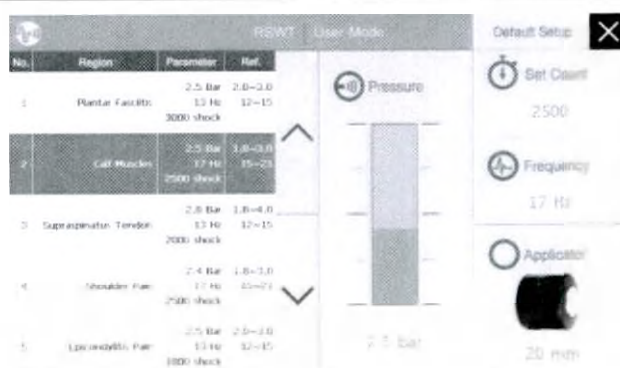
(f) The Parameter Settings screen appears, allowing the user to change parameter values.

(g) Press Region of the mode for 2~3 seconds and a screen will appear in which a name of the region can change. After changing the name, press Exit to close the screen.

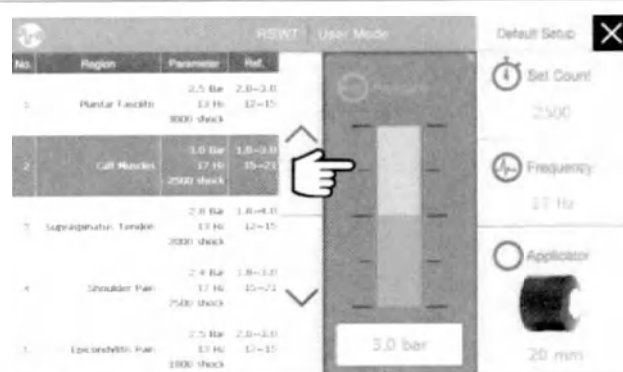(h) Set Count, Pressure, Frequency can change by selecting them and turning the encoder to change the values.

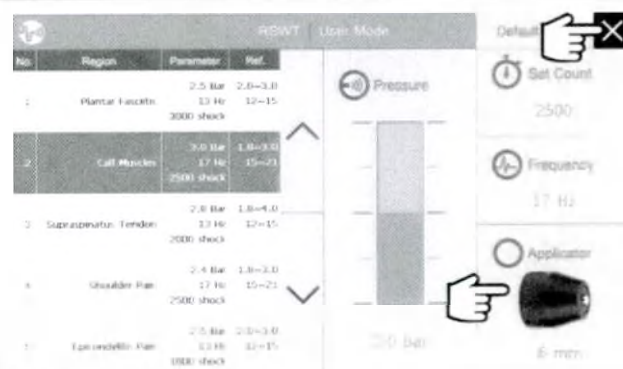(i) The user presses the applicator on the screen repeatedly until the type of the applicator required is displayed. When the user has finished setting the parameters, close the parameter setting screen by clicking close icon at the upper right corner.

(j) The changed parameter value and the selected applicator can be seen at the bottom of the mode window, in the parameter window, and in the applicator display of the channel selection.

(k) If you press the start button in the display once and then press the start button on the gun once, a shock wave will be generated (1 Shock). Press the start button on the gun for 2 to 3 seconds to start operation of the equipment. Pressure (Bar) and Frequency (Hz) can be changed during operation.

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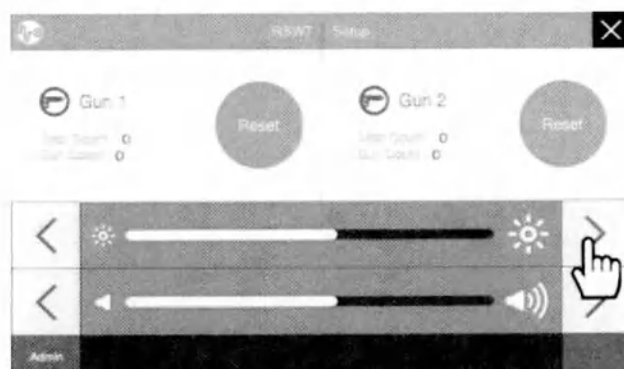
7. Messages

7.1 Over Temperature

If the air compressor in the equipment is overheated during operation, a warning as shown in [Fig. 13] will appear and the generation of shockwave will be temporarily stopped. In this case, do not turn off the power of the equipment, wait for a while until the message disappears while the power is on, and it will be re-enabled. If you see this message frequently, please contact the REMED Customer Service Center for a checkup.

Figure 13. Over temperature message

7.2 Gun Open

This equipment always checks the connection status of the gun. When the gun corresponding to the channel in use is not connected, the message "Gun Open" appears on the main screen as shown in Figure 14 below.

Figure 14. Gun open message

7.3 Components exchange messages

Depending on the number of shocks that are output, this equipment shows you the message for replacement of components as shown in Figure 15, 16. Even if the message appears, the equipment works, but it is recommended to replace the components for safe use of the equipment.

(1) Tube & Projectile exchange

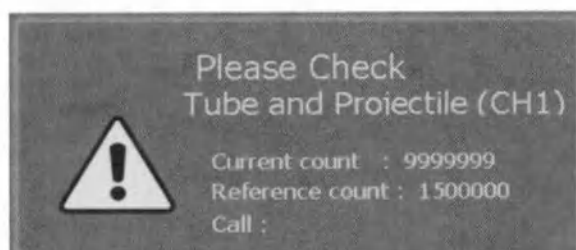
Figure 15. Tube & Projectile exchange message

(2) Handpiece (Gun 1/2) exchange

Figure 16. Handpiece (Gun) exchange message

8. Treatment protocols

The available treatment protocols of the SALUS-RSWT are shown in the Table below. Please note that this mode may vary by individual and pain level for reference purposes

Table 1. Treatment protocols

Application	Pressure (Bar)	Mode
Plantaritis	2.5	1
Pelvic pain	2.5	2
Superficial tendinitis	2.8	3
Shoulder pain	2.4	4
Perineal pain	2.4 ~ 2.5	5, 6
Achilles tendon pain	2.4 ~ 2.5	7, 8
Patellar tendon	2.0	9, 10
Tibial pericarditis	2.4	11, 12
bursitis	1.8 ~ 3.2	13, 14
Mucous membrane	2.0	15, 16
Cubital	3.0	18
Thigh pain	2.0	19
Stomach pain	2.0	20

- ☐ On the day of treatment, a patient may feel that the symptoms are getting worse, but the next day the patient may feel better.
- ☐ The intensity of treatment may vary slightly depending on the patient. The user must adjust the treatment intensity to suit patients.

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The SALUS-RSWT System is intended for use in the electromagnetic environment specified below. The customer or the user of the equipment should assure that it is used in such an environment.

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Samples that require EMC compliance: power cord, shock wave gun with cable, display, holder arm with transducer holder assembly, electromagnetic transducer with hose. ! The use of transducers and cables that differ from those installed is not allowed. This can result in increased electromagnetic emissions or reduced noise immunity. The OPTIMUS-Pro system should not be used in close proximity to or in conjunction with other equipment.

Emission Test	Compliance	Electromagnetic environment - guidance
RF emissions EN55011	Group 1	The Salus- RSWT system is intentionally generated and/or used conductively coupled radio-frequency energy which is necessary for the internal functioning of the equipment itself.
RF emissions EN55011	Class A	The Salus- RSWT system is suitable for use in all establishments other than domestic and those directly connected to the public low voltage power supply network that supplies buildings used for domestic purposes.
Harmonic Emissions IEC 61000-3-2	Class A	
Voltage fluctuations/ flicker emissions IEC 61000-3-3	Complies	

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The Salus- RSWT System is intended for use in the electromagnetic environment specified below. The customer or the user of the equipment should assure that it is used in such an environment.

Immunity Test IEC 60601	IEC 60601 Test Level	Compliance Level	Electromagnetic environment - guidance
Electrostatic discharge (ESD) IEC 61000-4-2	±6 contact ±8 air	±6 contact ±8 air	Floors should be wood, concrete or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30%.
Electrical fast Transient/ burst IEC 61000-4-4	±2 kV for power supply lines	±2 kV for power supply lines	Mains power quality should be that of a typical commercial or hospital environment.
Surge IEC 61000-4-5	± 1 kV differential mode ± 2 kV common mode	± 1 kV differential mode ± 2 kV common mode	Anti-surge protection needs to be incorporated into the main supply to the equipment if surge protection is to be guaranteed.
Voltage dips, short interruptions and voltage variations on power supply input lines. IEC 61000-4-11	<5 % U_T (>95 % dip in U_T) for 0,5 cycle 40 % U_T (60 % dip in U_T) for 5 cycles 70 % U_T (30 % dip in U_T) for 25 cycles <5 % U_T (>95 % dip in U_T)	<5 % U_T (>95 % dip in U_T) for 0,5 cycle 40 % U_T (60 % dip in U_T) for 5 cycles 70 % U_T (30 % dip in U_T) for 25 cycles <5 % U_T (>95 % dip in U_T)	Mains power quality should be that of a typical commercial or hospital environment. If the user of the SALUS-RSWT System requires continued operation during power interruptions, it is recommended that the SALUS-RSWT System be powered from an interruptible power supply.

	for 5 sec	for 5 sec	
Power Frequency (50/60Hz) Magnetic Field IEC 61000-4-8	3 A/m	3 A/m	Power frequency magnetic fields should be at levels characteristic of a typical location in a typical commercial or hospital environment
Radiated RF Immunity EN61000-4-3	80 MHz – 2.5 GHz 2Hz 80% amplitude modulation		Equipment should only be used in the vicinity of other equipment compliant with EN60601-1-2.
Conducted RF Immunity EN61000-4-6	0.15 MHz – 80 MHz 2Hz 80% amplitude modulation		Equipment should only be used in the vicinity of other equipment compliant with EN60601-1-2.

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The SALUS-RSWT System is intended for use in the electromagnetic environment specified below. The customer or the user of the equipment should assure that it is used in such an environment.

Immunity Test	IEC 60601 Test Level	Compliance Level	Electromagnetic environment - guidance
Conducted RF IEC 61000-4-6	3 V _{rms} 150 kHz – 80	3 V _{rms}	Portable and mobile RF communications equipment should be used no closer to any part of the SALUS-RSWT System, including cables, than the recommended separation distance calculated from the equation applicable to the frequency transmitter.
Radiated RF IEC 61000-4-3	3 V/m 80 MHz – 2.5	3 V/m	
			Recommended separation distance $d = \left[\frac{3.5}{V_1} \right] \sqrt{P}$ $d = \left[\frac{3.5}{E_1} \right] \sqrt{P}$ 80 MHz – 800 MHz $d = \left[\frac{7}{E_1} \right] \sqrt{P}$ 800 MHz – 2.5 GHz Where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer and d is the recommended separation distance in meters (m). Field Strengths from fixed RF transmitters, as determined by an electromagnetic site survey^a should be less than the compliance level in each frequency range^b.

		Interference may occur in the vicinity of equipment marked with the following symbol:
--	--	--

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a. Field strengths from fixed transmitters, such as base stations for radio (cellular/ cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the location in which the SALUS-RSWT System is used exceeds the applicable RF compliance level above, the SALUS-RSWT System should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as re-orienting or relocating the SALUS-RSWT System

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The SALUS-RSWT System is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. The customer or the user of the SALUS-RSWT System can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters) and the SALUS-RSWT System as recommended below, according to the maximum output power of the communications

Rated maximum output power of transmitter, W	Separation distance according to frequency of transmitter		
	150 kHz - 80 MHz $d = \left[\frac{3.5}{V_1} \right] \sqrt{P}$ Where $V_1 = 3$	80 MHz - 800 MHz $d = \left[\frac{3.5}{E_1} \right] \sqrt{P}$ Where $E_1 = 3$	800 MHz - 2.5 GHz $d = \left[\frac{7}{E_1} \right] \sqrt{P}$ Where $E_1 = 3$
0.01	0.117	0.117	0.233
0.1	0.369	0.369	0.738
1	1.167	1.167	2.333
10	3.689	3.689	7.379
100	11.667	11.667	23.333

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NOTE 2. These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people..

11. FAQ

Q1) Is the applicator semi-permanent?

- A) Applicator should be purchased as a consumable item. Please be careful because it can shorten the life time when you use it without contact with the affected part of the body.
- B) Do not use it with very high intensity (pressure) because it may shorten the life time. Please use it with appropriate intensity.

Q2) When should the projectile be replaced?

- A) It guarantees over 1,500,000 with good maintenance. It is recommended to clean it regularly.

Q3) Do I have to use a gel?

- A) It is recommended to use it to deliver the shockwave without loss to the affected area of the body.

Q4) How do I maintain the equipment?

- A) Projectile guide (pipe) cleaning is basically recommended to clean the oil at least once a week. If you use it a lot, we recommend cleaning the same day. Please refer to 12.3 Cleaning procedure.

Q5) Can I use 2 guns (channels) simultaneously?

- A) We provide one gun with the equipment and it is not possible to use two guns simultaneously.

12. Maintenance

12.1 Routine inspection of equipment

- ☐ The covering of power line of equipment, gun connecting line, etc. shall not be peeled off and internal lines shall not be exposed, and shall not be damaged by impact from outside.
- ☐ Wash the outside of equipment so that there is no foreign material.
- ☐ The button for equipment operation, etc. must not shake.
- ☐ The various parts attached to equipment must not shake.
- ☐ If any of the above occur, contact customer service for help

12.2 Safety inspection

- ☐ In order to ensure safe use, internal cleaning should be performed once per year by a person authorized by this company.
- ☐ In order to ensure safe use, be sure to check the equipment including internal components and output voltage from the person who has been given authorization from the company once per year.
- ☐ Please clean the gun before storing it
- ☐ When storing the product for a long period of time, be sure to check the product before using it.
- ☐ Please note the followings regarding storage conditions.
 - Keep out of water
 - Keep away from direct sunlight
 - Do not store near heaters
 - Avoid locations subject to excessive shock or vibration, exposure to chemicals or explosive gases.

12.3 Cleaning procedure

12.3.1 General cleaning

The SALUS-RSWT and its accessories can be cleaned in variable ways. Avoid any damage or contamination of the equipment using the methods recommended below. Before cleaning the unit, be sure to turn off the unit. Cleaning of the inside of the machine is only allowed to the person authorized by us.

- ☐ Periodically clean the equipment's exterior surface, gun and touch screen with a soft cloth moistened with alcohol. Do not use abrasives, lacquers, thinners, ethylene, or oxides, as these materials may cause permanent damage to the equipment.
- ☐ Do not submerge any section of equipment into liquid or detergent. In addition, no liquid must be allowed to enter the equipment or accessories.
- ☐ When the operation for one patient is completed, clean the applicator with alcohol and smooth cloth.
- ☐ If you want to use a disinfectant or a method not mentioned above, it may damage the device. Please note that in the event of damage to the equipment using unauthorized materials, no warranty of any kind is warranted.

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For cleaning, use running water to clean. If it is hardly fixed, it is effective to immerse the applicator in water for 2 to 3 hours and then wash it. After washing, remove the water using a soft dry cloth and dry it at room temperature for over 1 hour. Remove the applicator as shown in Figure 17 and clean the front metal part of the gun with a soft brush

Figure 17. Applicator cleaning

(1) 6 mm Applicator cleaning

Disassemble the applicator from the gun as shown in Figure 18 below and clean the gel adhered to the arrow in the below.

Figure 18. 6 mm Applicator cleaning

(2) 10 mm Applicator cleaning

Disassemble the applicator from the gun as shown in Figure 19 and clean the gel adhered to the arrow in the below.

Figure 19. 10 mm Applicator cleaning

(3) 15 mm Applicator cleaning

Disassemble the applicator from the gun as shown in Figure 20 below and clean the gel adhered to the arrow in the below.

Figure 20. 15 mm Applicator cleaning

(4) 20 mm Applicator cleaning

Disassemble the applicator from the gun as shown in Figure 21 and clean the gel adhered to the arrow in the below.

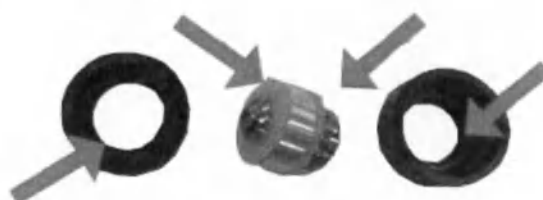
Figure 21. 20 mm Applicator cleaning

(5) 35 mm Applicator cleaning

Disassemble the applicator from the gun as shown in Figure 22 below and clean the gel adhered to the arrow in the below.

Figure 22. 35 mm Applicator cleaning

12.3.3 Projectile & projectile guide cleaning

As shown in Figure 23-(1) rotate the applicator to disassemble and rotate the center part of the gun to disassemble. Use the tweezers to remove the projectile guide in the center of the separated gun or the projectile left inside the gun, Figure 23-(2). Clean the projectile using a dry rag, and clean the inside of the gun using the cleaning tools provided Figure 23-(3).

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Figure 23. Projectile & projectile guide cleaning

12.3.4 Connection part of the gun cleaning

After a long period of use, a small amount of water and oil may flow out from the areas indicated by the arrows as shown in Figure 24. This is not a malfunction, so you can wipe it off with a soft cloth. However, if this happens frequently, please contact the REMED Customer Service Center for inspection.

Figure 24. Gun connection part cleaning

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		
Projectile	Projectile guide	Bumper

Figure 25. Projectile, projectile guide, bumper

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(2) Exchanging method

As shown in Figure 26-(1) rotate the applicator to disassemble and rotate the center of the gun to disassemble. Use the tweezers to remove the projectile guide in the center of the separated gun or the projectile left inside the gun, Figure 26-(2). The bumper can be also removed by using the tweezers. Use hex wrench to remove the projectile guide as shown in Figure 26-(3). After replacing new parts, assemble in reverse order.

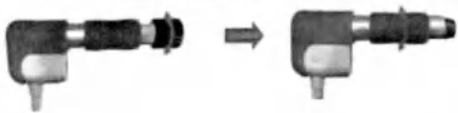
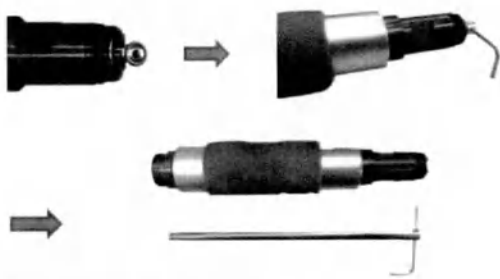
(1) Disassemble the applicator and gun.	
(2) Use the tweezers to remove the projectile guide in the center of the separated gun or the projectile left inside the gun.	
(3) Remove the bumper using the tweezers.	
(4) Remove the projectile guide using the hex wrench.	

Figure 26. Projectile, projectile guide, bumper disassembling

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(2) Method of condensate remove

There is a black butterfly valve below the front of the equipment as shown in Figure. 27]. Turn on the power and turn the valve in the right (clockwise) direction while the pressure in the tank is raised

Figure 27. Air compressor condensate remove

When the condensate is removed, opening the valve can cause oil to come out and contaminate the environment. Combine the transparent hose provided with the valve and prepare the container to receive condensate. Then open the valve to remove condensate. When connecting the transparent hose to the valve, insert the hose into the blue one-touch fitting sleeve, Figure 28.

Please be careful because it may splash in all directions when condensate is removed. After all of the condensate is coming out, remove the hose by pressing the sleeve of the one-touch fitting as shown in Figure 29.

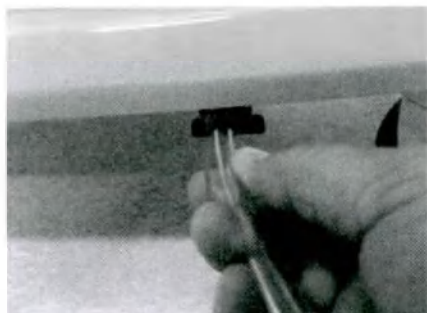
Figure 28. Hose connection

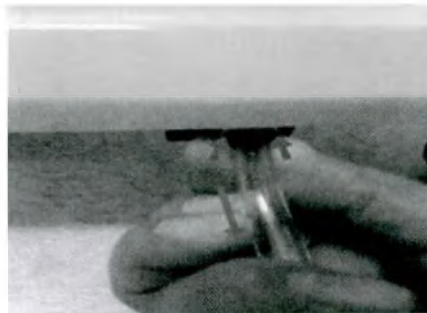
Figure 29. Hose disconnection

**Caution**

After removing the condensate, remove the transparent hose from the one-touch fitting, and then turn the butterfly valve to the left (counter-clockwise) direction to lock it.

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Table 2. Self-troubleshooting

Symptom	What to do	References in the Instructions For Use
Equipment does not turn on.	Check if the power connector of the equipment is properly connected.	• 3.4.2 Figure 3. Power cable connection
	Check if the power switch of equipment is turned on.	• 4.1.2 Figure 8. Rear view
	Confirm if the Key Lock switch ON.	• 4.1.1 Figure 7. Front view
	Confirm if the emergency button on the equipment is pressed. If it is pressed, turn it to the right to release it	• 4.1.1 Figure 7. Front view
Over Temperature message is displayed.	Air compressor overheat message. Wait a moment and the message will disappear and be re-usable. If you receive this message frequently, please	• 7.1 Over temperature

Symptom	What to do	References in the Instructions For Use
	contact the REMED Customer Center for a checkup.	
	If the room temperature is too high, disorder can be caused in the cooling. Maintain the room temperature at less than 25°C as much as possible.	• 3.3 Environment condition
When the gun is in operation, there is a windy sound.	Check the air tank valve lock. Turn the valve counterclockwise to lock it.	• 12.5 Figure 27. Air compressor condensate remove
	Please check the connection of gun connection cable..	• 4.1.2 Figure 6. Rear view
Gun's output is weak.	Moisture in the air tank can cause irritation. Please check and remove condensate.	• 12.5 Removing air compressor condensate
	If there is a dust inside the gun, or the applicator is loosened slightly, the output may be weak.	• 12.3.3 Projectile & projectile guide cleaning
There is no output from the gun.	Please check the connection of gun connection cable..	• 4.1.2 Figure 6. Rear view
	After you press Start on the touch screen, you must press the Start on the gun to get the output.	• 6.2 Operation mode—(k)
	It may be relating the maintenance. (1) Replace the projectile guide and projectile. (2) Remove condensate (3) Check applicator assembly status Please.	• 12.4. Replacement of components • 12.5 Removing air compressor condensate
Wheel does not roll when moving the equipment.	Move the equipment by placing the wheel in the "Loose" position..	• 3.4.3 Figure 4. Transporting wheel in "Loose" and 'Locked'.

Symptom	What to do	References in the Instructions For Use
In the following cases, stop the operation by cutting off the power to the equipment, and contact the service center. → The main power switch spontaneously turns off. → The LCD screen of operation panel does not illuminate when power is turned off and then turned on again. → When stimulation is not generated by gun. → The temperature icon on the screen blinks and the equipment is not operated. → The message indicating no connection between the gun and the main unit appears constantly.		

12.7 Regular inspection for performance maintenance

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12.8 Warranty

- ☐ This product is manufactured through an intensive quality control and inspection process.
- ☐ Compensation for repair and exchange of the product will be in accordance with the "Consumer Injury Compensation Rule" announced by the Economy Planning Board.
- ☐ The warranty period for this equipment is 1 year.
- ☐ If a failure occurs during the warranty period under normal operating conditions, the equipment will be repaired free of charge by our service center.
- ☐ If a problem with the equipment occurs during the warranty period, prepare the following customer service request to report to our company.

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- Name of model : SALUS-RSWT
- Date of purchase :
- Serial number of product :

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appeared before me and admitted said principal's
subscription to the attached

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User Manual.
This is hereby attested on this
26th day of Jul. 2023 at this office.

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Daehan Law Firm & Notary Office Inc.
Belong to Seoul Central District
Prosecutors' Office
DaeGwang B/D 3F, 36 Dongjak-daero,
Seocho-gu, Seoul, Korea

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Attorney-in-law

Handwritten signature of Kwon Hyuk Jung

KWON HYUK JUNG

This office has been authorized by the
Minister of Justice, the Republic of
Korea, to act as Notary Public since
January 6, 2010 under Law No. 11823

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
TEL : (02) 590-0900(代)

FAX : (02) 3477-0957

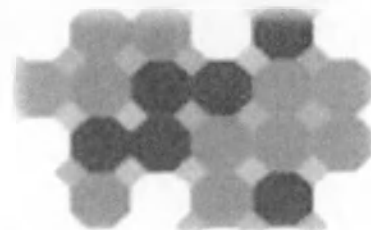REMED
Rehabilitation Equipment Company

User Manual

Device for extracorporeal shock wave therapy (ESWT),
model SALUS -OPTIMUS-Pro

Remed Co., Ltd.

Co-CEO
Eun-Hyun, Koh Sang-Yong, Lee

REMED Co., Ltd., Republic of Korea

Preface

The user of this device must sufficiently understand its functions and the precautions that must be observed for its safe use and stable performance.

For the safe use and post management of the device, be sure to thoroughly understand the details of the Instructions For Use before using the device. This is essential to guaranteeing the safe use and stable performance of the device.

The Instructions For Use provides a guide to the efficient use of the OPTIMUS-Pro. For clinical definitions, pathological effects etc., of its functions, refer to the related medical publications.

As the Instructions For Use has been arranged in independent chapters, some of the descriptions have been duplicated.

If you encounter any problems during operation of the equipment, stop its use immediately, and contact the customer service center of REMED (Refer to chapter 12).

Copyright

This Instructions For Use is protected by the Copyright REMED.

Transferring or storing all or part of the contents of the Instructions For Use to another media through translation, copying, scanning, photography, etc. is prohibited without the prior written approval of REMED.

The Instructions For Use may include editing errors. In addition, technical changes to improve performance may be made to this device after this manual has been published, without any prior notice.

Contents

1. Information on safety	10
1.1 Definitions of symbols	10
1.2 Requirements for safety	12
1.1.1 Precautions for electrical safety	12
1.1.2 Classifications	12
1.3 Precautions for Use	13
1.4 Contraindication	15
1.4.1 Extracorporeal Shock Wave Therapy	15
1.4.2 Electromagnetic Stimulator	15
1.5 Side effect	17
1.5.1 Extracorporeal Shock Wave Therapy	17
1.5.2 Electromagnetic Stimulator	17
1.6 General warnings	17
1.7 Safety Label	17
2. Intended use	18
2.1 Extracorporeal Shock Wave Therapy	18
2.2 Electromagnetic Stimulator	18
3. Installation	19
3.1 Components	19
3.2 Electrical installation condition	20
3.3 Environmental condition	20
3.3.1 Operation environment	20
3.3.2 Transport and Storage environment	20
3.4 Installation method	21
3.4.1 Precautions during installation	21
3.4.2 Connection of power	21
3.4.3 Moving and fixing of equipment	22
3.4.4 Transducer	20
3.4.5 Gun 21	
4. Product description	22
4.1 External view	22

4.1.1 Front view	22
4.1.2 Rear view	22
4.2 Transducer and gun	23
4.3 Controls	23
4.3.1 Electromagnetic field function control	23
4.3.2 Extracorporeal shockwave function control	24
5. Preparations for operation	25
5.1 Moving and fixing equipment	25
5.2 Power connection	25
6. Operation instructions	26
6.1 Transducer placement	26
6.2 Applications	27
6.3 Operation mode	29
6.3.1 Magnetic stimulation function	29
6.3.2 Extracorporeal shockwave function	35
6.3.3 Environment setting	39
7. Messages	40
7.1 Over Temperature - Transducer	40
7.2 Transducer Error	40
7.3 Over Temperature - Air compressor (RSWT)	41
7.4 Gun Open	41
7.5 Components exchange messages	42
8. Technical specifications	43
8.1 Size and weight	43
8.2 Power	43
8.3 Modes – Magnetic stimulation	43
8.4 Modes - RSWT	43
8.5 RSWT number of shocks	43
8.6 RSWT intensity	43
8.7 RSWT frequency	43
8.8 RSWT applicator	43
8.9 Magnetic field intensity	43

8.10	Life Time.....	44
9.	EMC emissions and immunity	50
10.	FAQ	56
11.	Maintenance	57
11.1	Routine inspection of equipment	57
11.2	Safety inspection	57
11.3	Cleaning procedure	58
11.3.1	General cleaning.....	58
11.3.2	Cleaning applicators	59
11.3.3	Projectile & projectile guide cleaning	61
11.3.4	Connection part of the gun cleaning	61
11.4	Replacement of components	60
11.5	Air compressor condensate remove	61
11.6	Troubleshooting.....	62
11.7	Regular inspection for performance maintenance	64
11.8	Warranty	64
12.	Contact information.....	66

Figure Contents



Figure 1. Power cable connection	21
Figure 2. Transporting wheel in 'Loose' (left) and in 'Locked' (right)	22
Figure 3. Separating the transducer from the transducer arm supporter	20
Figure 4. Adjusting height of the transducer arm supporter	21
Figure 5. Front view	22
Figure 6. Rear view	22
Figure 7. Transducer and gun	23
Figure 8. Controls - Magnetic stimulation	23
Figure 9. Controls - RSWT	24
Figure 10. Changing the angle of the transducer	26
Figure 11. Deflected transducer	26
Figure 12. Application of the magnetic stimulation	27
Figure 13. Application of the RSWT	28
Figure 14. REMED logo Figure 15. Function selection screen	29
Figure 16. Magnetic stimulation main screen	29
Figure 17. RSWT main screen	35
Figure 18. Over temperature messages – Transducer	40
Figure 19. Transducer error message	41
Figure 20. Over temperature - Air compressor (RSWT)	41
Figure 21. Gun open message in RSWT function screen	42
Figure 22. Gun open message in magnetic stimulation function screen	42
Figure 23. Tube & Projectile exchange message	42
Figure 24. Gun exchange message	42
Figure 25. Applicator cleaning	59
Figure 26. 6mm Applicator cleaning	59
Figure 27. 10mm Applicator cleaning	59
Figure 28. 15mm Applicator cleaning	59
Figure 29. 20mm Applicator cleaning	60
Figure 30. 35mm Applicator cleaning	60
Figure 31. Projectile & projectile guide cleaning	61
Figure 32. Gun connection part cleaning	61
Figure 33. Projectile, projectile guide, bumper	60
Figure 34. Projectile, projectile guide, bumper disassembling	60
Figure 35. Air compressor condensate remove	61
Figure 36. Hose connection Figure 37. Hose disconnection	61

.. Information on safety

1.1 Definitions of symbols

- When there is specific information in the Instructions For Use that needs to be emphasized for safety, the following terms and symbols are indicated. All warnings and precautions shall always be observed.
- The manufacturer or agent of the product is not responsible for any personal/material damage caused by erroneous use, operation for purposes other than its intended objective and negligence of product maintenance.

	Warning	"Warning" symbol is used to indicate a life-threatening risk to the operator if the warning is disregarded.
	Caution	"Caution" symbol is used to indicate that injury or damage can be caused if the caution is disregarded.
	No pushing	To prohibit pushing against an object
	No sitting	To prohibit sitting on a surface
	BF Type	Applied part BF type)
	Operation instructions	Operation instructions
	instructions for use.	Follow instructions for use.
	Serial number	Serial number
	Article number	Article number
	Manufacturer	Manufacturer
	Date of manufacture	Date of manufacture
	Do not stack	To indicate that the items shall not be vertically stacked, either because of the nature of the transport packaging or because of the nature of the items themselves.

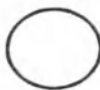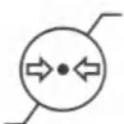	Keep away from rain	To indicate that the transport package shall be kept away from rain and in dry conditions.
	Fragile	To indicate that the contents of the transport package are fragile and the package shall be handled with care.
	This way up	To indicate correct upright position of the transport package.
	Recycling	To indicate the location of a recycling bin or container.
	"OFF" (Power)	To indicate disconnection from the mains, at least for mains switches or their positions, and all those cases where safety is involved.
	"ON" (Power)	To indicate connection to the mains, at least for mains switches or their positions, and all those cases where safety is involved.
	Alternating current	To indicate on the rating plate that the equipment is suitable for alternating current only to identify relevant terminals.
	WEEE	Indicates that when the end-user wishes to discard this product, it must be sent to separate collection facilities for recovery and recycling.
	Temperature limitation	To indicate the maximum and minimum temperature limits at which the item shall be stored, transported or used.
	Humidity Limitation	To indicate the acceptable upper and lower limits of relative humidity for transport and storage.
	Atmospheric pressure limitation	To indicate the acceptable upper and lower limits of atmospheric pressure for transport and storage.
	Dangerous Voltage	Warning, electricity

1.2 Requirements for safety

1.1.1 Precautions for electrical safety

- ❑ The rated power for this equipment is 220 V~.
- ❑ Confirm that all the connecting sections (power line or peripheral equipment) are adequately connected to the equipment.
- ❑ Confirm that the equipment is completely grounded.
- ❑ Repair, expansion, and installation of equipment shall not be performed by anyone other than the specialized personnel authorized by the manufacturer. Arbitrary disassembling/assembling of equipment by the user is absolutely prohibited.
- ❑ Before connecting other equipment that is not specified in the Instructions For Use, be sure to inform this company or the agency with authority over product marketing.
- ❑ To avoid electrical noise during use, the equipment shall be installed at a significant distance from any generator, X-ray equipment, broadcasting device, mobile electric wire, and so forth.
- ❑ Needs attention due to potential electromagnetic or other interference between ME EQUIPMENT and other devices.
- ❑ An independent power circuit is essentially required, and sharing of power with other electronic devices is not recommended.
- ❑ Do not place the ME EQUIPMENT in a location where it is difficult to disconnect the power plug or other detachable plug.

1.1.2 Classifications

- ❑ Protection type and level for electric shock: Class I, BF Type
- ❑ Electromagnetic compatibility (EMC) test standard: class A

1.3 Precautions for Use

- ❑ Installation and re-installation of equipment shall always be performed by specialized personnel authorized by the manufacturer.

	Caution	"PLEASE MAKE SURE THAT PATIENT IS NOT WEARING ANYTHING WITH METALLIC OBJECTS DURING TREATMENT FAILURE TO DO SO MAY RESULT IN SERIOUS SKIN BURN!"
---	----------------	--

- ❑ Installed equipment shall be subject to regular safety inspections from specialized personnel authorized by the manufacturer.
- ❑ Repair or installation of equipment shall only be performed by the specialized personnel authorized by the manufacturer, and as such, all arbitrary disassembly/assembly of equipment by the user is absolutely prohibited.
- ❑ For electrical safety, equipment shall always be connected to a safely grounded power supply for operation.
- ❑ The operator and manager of equipment shall thoroughly understand the Instructions For Use, and it must be kept in the close proximity of the equipment.
- ❑ To prevent safety accidents and ensure adequate maintenance, place the safety marks, guide phrases and regular inspection table provided with the equipment at an easily visible location near the equipment.
- ❑ As a fan is installed in the body of the equipment for air circulation, remove curtains or any other objects which can block the flow of air from near the equipment.
- ❑ From the surroundings of the location where equipment is installed, remove water, alcohol, flammable materials, and so forth.
- ❑ As a strong magnetic field is generated around the magnetic field generating section, equipment operation technicians, assistants, and patients must not hold any items which can be affected by the magnetic field.

	Note	Items such as wristwatches, mobile phones may be damaged by magnetic fields, so please be careful and keep them separately.
---	-------------	---

- ❑ When equipment is operated, do not use any mobile phones, radio sets, mobile radio transmitters, wireless toys, etc. in the proximity of the equipment.
- ❑ During the operation of the equipment, the patient shall not take drinks, water, etc. which can influence the equipment.
- ❑ Do not separate gun and applicator during operation of the equipment, as it may cause injury due to recoil of airpressure.
- ❑ If you do not apply the Transmission Gel enough or do not come into contact with the skin, patients may feel pain because the air bubble is formed and the shock wave is

generated. Therefore, it is necessary to apply gel sufficiently and close the gun to patients tightly to prevent pain.

- ❑ Shockwave may be affected by overpressure when used around the lungs, ear, and gastrointestinal tract.
- ❑ Do not directly apply shockwaves to the heart, head or wounded area.
- ❑ This equipment should be used for patients who can accurately express and communicate pain during operation.
- ❑ Excessive pressure on the stimulation area with a gun may cause skin damage such as blisters and redness.
- ❑ Do not move the unit during use because the gun may fall or the unit may be damaged due to impact or movement of the unit.

1.4 Contraindication

1.4.1 Extracorporeal Shock Wave Therapy

Do not use in combination with other electronic medical devices.
Do not use to carry out procedures in the area of the endoprosthesis.

The following patients must be approved by the physician during the procedure:

- Patients with high body temperature (more than 37 °C), the elderly (over 85 years old) and children (under 3 years old)
- a person with implants (for example, an artificial pacemaker, an implantable infusion port, a hearing aid).

Do not treat with the product the following patients who have

- blood disorders, clotting problems or use of anticoagulants;
- blood-thinning drugs (warfarin);
- thrombosis;
- oncological diseases;
- with diabetes mellitus, diabetic neuropathy;
- pregnant women;
- acute inflammation in the target area;
- growth of cartilaginous tissue in children;
- contact infectious patients (to avoid transmission of infection);
- Corticoid therapy for 6 weeks.

※ Caution: shock waves should not be applied on top of air-filled tissues (lungs), nor should they be applied near large nerves, vessels, vertebrae or heads.

1.4.2 Electromagnetic Stimulator

- ❑ Do not use the equipment in parallel with other electronic medical equipment.
- ❑ Be careful to ensure that magnetic stimulation does not penetrate the heart region.
- ❑ In general, patients in the following categories cannot be treated with this equipment. Prior to any treatment with this equipment, permission of the doctor in charge must be obtained.

- patients with high body temperature (more than 37 °C), pregnant women, the elderly (over 85 years old) and children (under 3 years old)
- a patient with a history or diagnosis of epilepsy or convulsive disability;
- a patient with suspected epilepsy based on electroencephalography;
- patients with traumatic brain injury and spinal cord injury within a month from the moment of injury;
- patients with implantable infusion port systems or hearing aids;
- patient with cranial implants.
- use during magnetic resonance imaging.

1.5 Side effect

1.5.1 Extracorporeal Shock Wave Therapy

- ☐ Swelling
- ☐ Redness
- ☐ Hematoma
- ☐ Viscous bleeding
- ☐ Pain

1.5.2 Electromagnetic Stimulator

- ☐ No side-effect

1.6 General warnings

- ☐ Be sure to use the equipment in accordance with the Instructions For Use.
- ☐ In the Instructions For Use, the method for the correct use of the equipment is described.
- ☐ Read the Instructions For Use carefully.
- ☐ The equipment shall be operated under the supervision of a person who has completed specialized medical education.
- ☐ This company is not responsible for any disadvantages / damages related to the operation of equipment by a person without a medical education.
- ☐ This equipment shall not be altered, remodeled, or used for any purpose other than the intended objective.

1.7 Safety Label

- ☐ Labels attached to devices provide information and safety information.

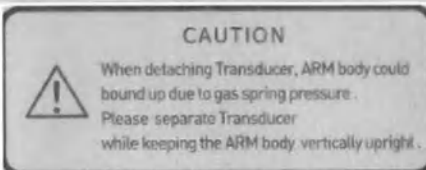
Label	Attachment position	Safety contents
 CAUTION When detaching Transducer, ARM body could bound up due to gas spring pressure. Please separate Transducer while keeping the ARM body vertically upright.	Upper side surface of transducer arm	Warning mark related to transducer arm

Table 1. Safety label

2. Intended use

The Optimus-Pro is a combo device of extracorporeal shockwave therapy and electromagnetic stimulator.

2.1 Extracorporeal Shock Wave Therapy

- ❑ The OPTIMUS-Pro is indicated for extracorporeal shock wave treatment of heel pain due to chronic proximal of heel pain due to chronic proximal plantar fasciitis for patients of age greater than 18 years with a history of failed alternative conservative therapies for at least 6 months. Chronic proximal plantar fasciitis is defined as traction degeneration of the plantar fascial band at the origin on the medial calcaneal tuberosity that has persisted for six months or more.

2.2 Electromagnetic Stimulator

- ❑ The OPTIMUS-Pro is an electrically powered device intended for medical purposes that repeatedly contracts muscles by passing electrical currents through electrodes non contacting the affected body area. And, in addition, the electro-magnetic stimulator is intended to provide entirely non-invasive electromagnetic stimulation of pelvic floor musculature for the purpose of rehabilitation of weak pelvic muscles and restoration of neuromuscular control for the treatment of urinary incontinence in human.

3. Installation

3.1 Components

Main body with transducer	Power cable	Instructions For Use										
		 <i>OPTIMUS-Pro Instructions For Use</i> <i>Extracorporeal Shock Wave Therapy & Electromagnetic Simulator</i>  <i>REMEDI Co., Ltd.</i> 										
Gun	Applicator											
	<table><tr><th>6mm</th><th>10mm</th><th>15mm</th><th>20mm</th><th>35mm(option)</th></tr><tr><td></td><td></td><td></td><td></td><td></td></tr></table>		6mm	10mm	15mm	20mm	35mm(option)					
6mm	10mm	15mm	20mm	35mm(option)								
												
Cleaning kit												
Cleaning stick  Tweezers  Hex wrench 												

Table 2. Main components

3.2 Electrical installation condition

- Input power : 220 V~, 50 Hz
- Power consumption : Max. 1.5 kVA

3.3 Environmental condition

3.3.1 Operation environment

- Temperature : 10 – 30 °C (50 – 86 °F)
- Humidity : 30 – 85 % RH
- Pressure : 70 kPa - 106 kPa

3.3.2 Transport and Storage environment

- Temperature : -10 – 60 °C (14 – 140 °F)
- Humidity : 80 % RH max
- Pressure : 70 kPa - 106 kPa

3.4 Installation method

3.4.1 Precautions during installation

- ❑ Install the equipment on a flat surface.
- ❑ Confirm that the power is connected.
- ❑ Operate the OPTIMUS-Pro within the ambient temperature range of 10 – 30 °C, and a humidity between 30 – 85 %.
- ❑ Install the equipment in a location with appropriate ambient temperature and humidity, and do not install in a location where it will be exposed to dust or flammable materials.
- ❑ Be careful not to damage the equipment by excessive shock.

3.4.2 Connection of power

- ❑ Confirm that power socket and power cable are grounded before operating the machine.
- ❑ Plug the power cable into the power terminal of the OPTIMUS-Pro as shown in below, and connect the other end to a grounded power socket.
- ❑ Confirm that the power cable is properly connected to the power terminal and power socket. If the connection is incorrect, unexpected problems can occur.
- ❑ Do not connect several devices to one power outlet.
- ❑ In particular, do not use a power cable that generates connection noise.

Figure 1. Power cable connection

Warning

To prevent the danger of electrical shock, connection shall be made to a protected and grounded power supply.

3.4.3 Moving and fixing of equipment

- ❑ Before moving the equipment, remove the connected power cable and the connection cables for peripheral equipment.
- ❑ When moving the equipment, place the transporting wheel at the lower section in the "loose" position.
- ❑ When the movement is completed, prevent the equipment from shaking by placing the transporting wheel in the "locked" position.
- ❑ If the equipment is moved while the wheel is in the "locked" position, it will damage the wheel, so please be careful.

Figure 2. Transporting wheel in 'Loose' (left) and in 'Locked' (right)

3.4.4 Transducer

- ❑ The OPTIMUS-Pro is equipment that generates a magnetic field by applying a strong current to a transducer. As heat is generated from the transducer due to the strong current used to create the magnetic field, cooling is performed by circulating the cooling oil inside the transducer. If the transducer is used while in a tilted position, the cooling oil will not reach all sections of the transducer, and injury or stoppage of operation due to overheating of the transducer can result. For this reason, it is recommended to place the transducer in a horizontal position as much as possible.

	Caution	Maintaining the transducer in a horizontal position is recommended.
---	----------------	---

- ❑ When separating the transducer from the transducer arm supporter, the pressure of the gas spring may cause the arm to jump up and the transducer may fall off. Be sure to hold the transducer with one hand while holding the arm supporter vertically. When replacing the transducer, be sure to mount the arm support vertically as well.

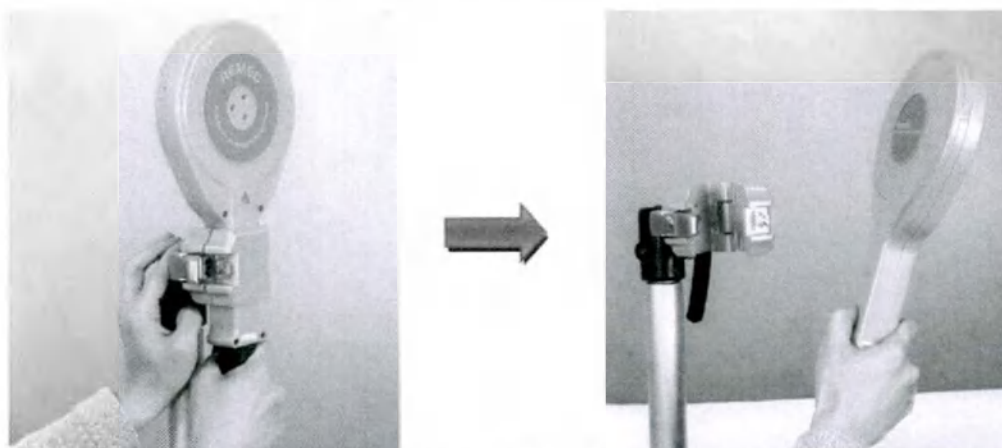
Figure 3. Separating the transducer from the transducer arm supporter

	Warning	When separating the transducer from the transducer arm supporter, the pressure of the gas spring may cause the arm to jump up and the transducer may fall off. Be sure to hold the transducer with one hand while holding the arm supporter vertically. When replacing the transducer, be sure to mount the arm support vertically as well.
---	----------------	---

- ❑ When using the OPTIMUS-Pro for a long time, the transducer arm support may be struck down or up. In this case, the position of the transducer arm support can be changed by adjusting the gas spring pressure. To lower the arm support, turn the

hexagonal wrench bolt in the (-) direction as shown in Figure 5, and turn it in the (+) direction to lift it upwards.

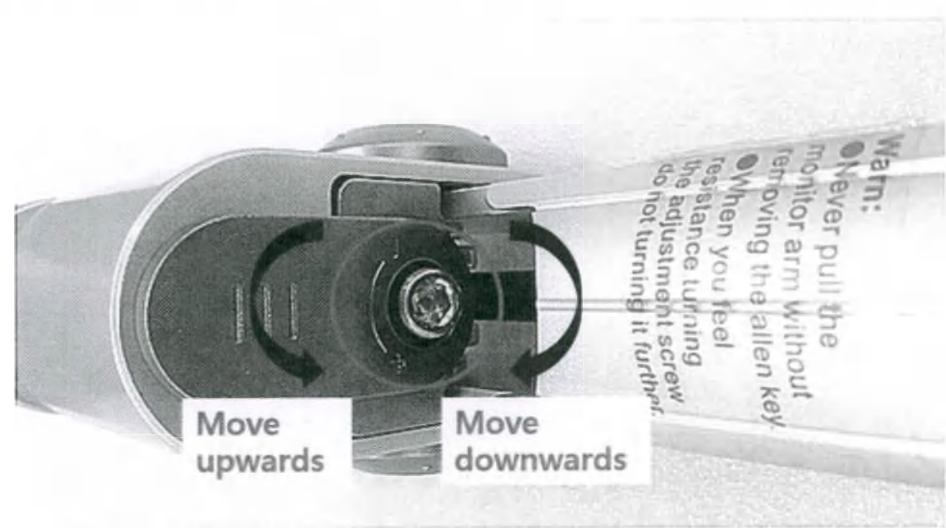
Figure 4. Adjusting height of the transducer arm supporter

3.4.5 Gun

Do not separate gun and applicator during operation of the equipment, as it may cause injury due to recoil of air pressure.

	Warning	Do not separate gun and applicator during the operation.
--	----------------	--

4. Product description

4.1 External view

4.1.1 Front view

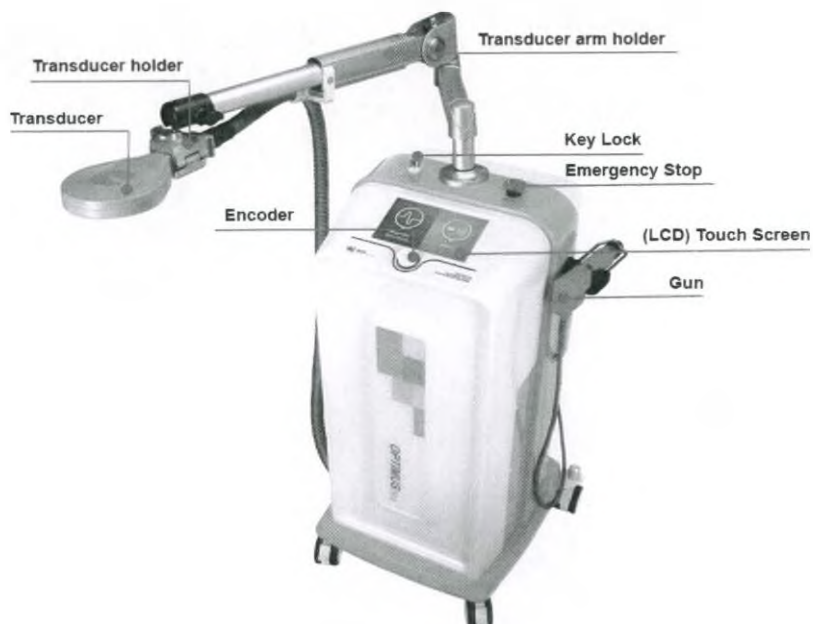
Figure 5. Front view

4.1.2 Rear view

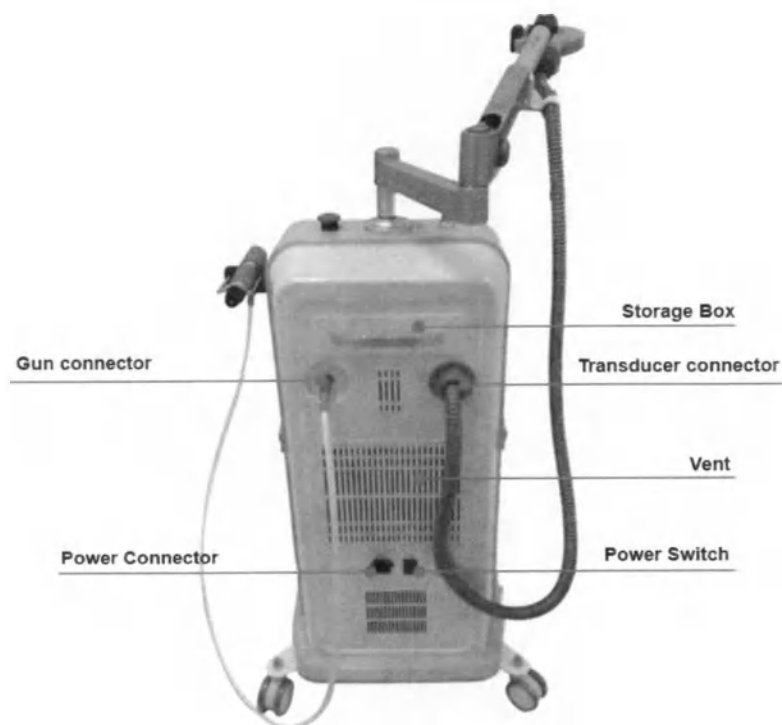
Figure 6. Rear view

4.2 Transducer and gun

Figure 7. Transducer and gun

4.3 Controls

4.3.1 Electromagnetic field function control

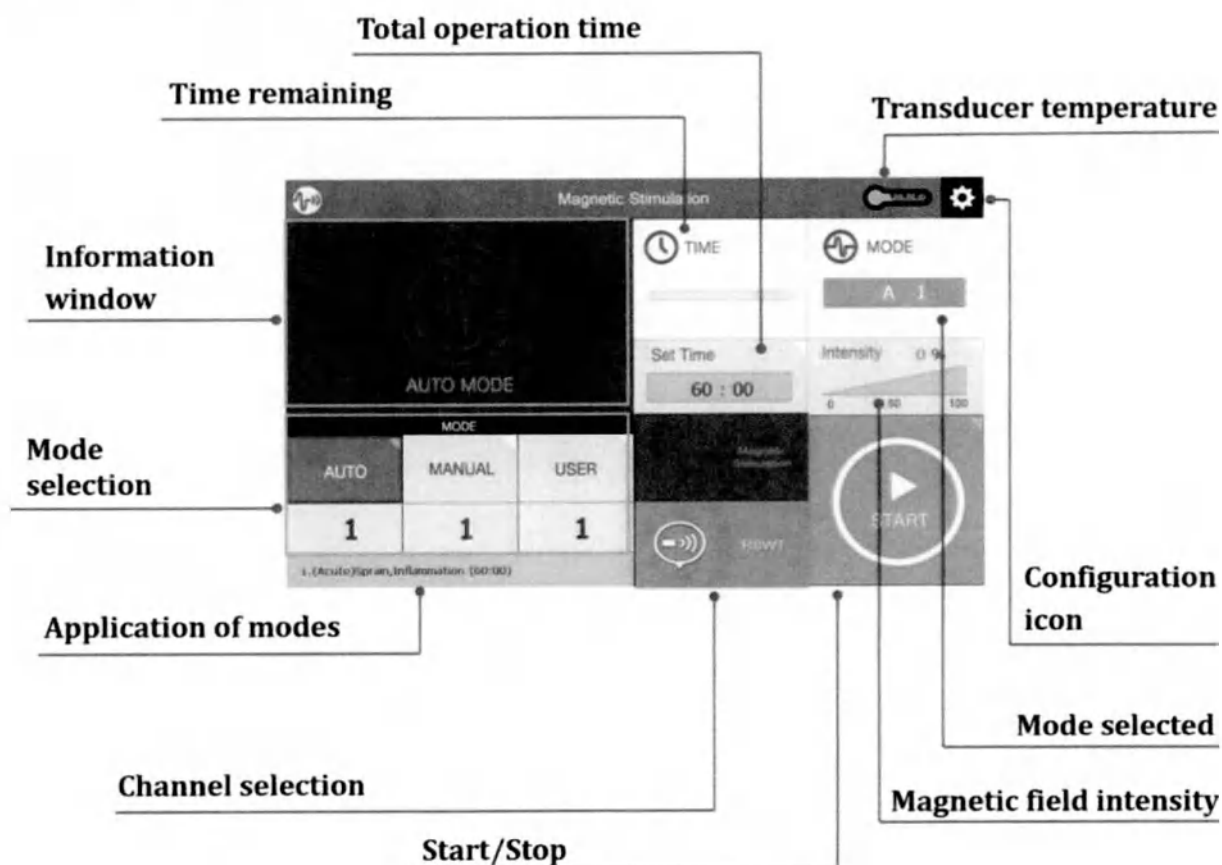
Figure 8. Controls - Magnetic stimulation

4.3.2 Extracorporeal shockwave function control

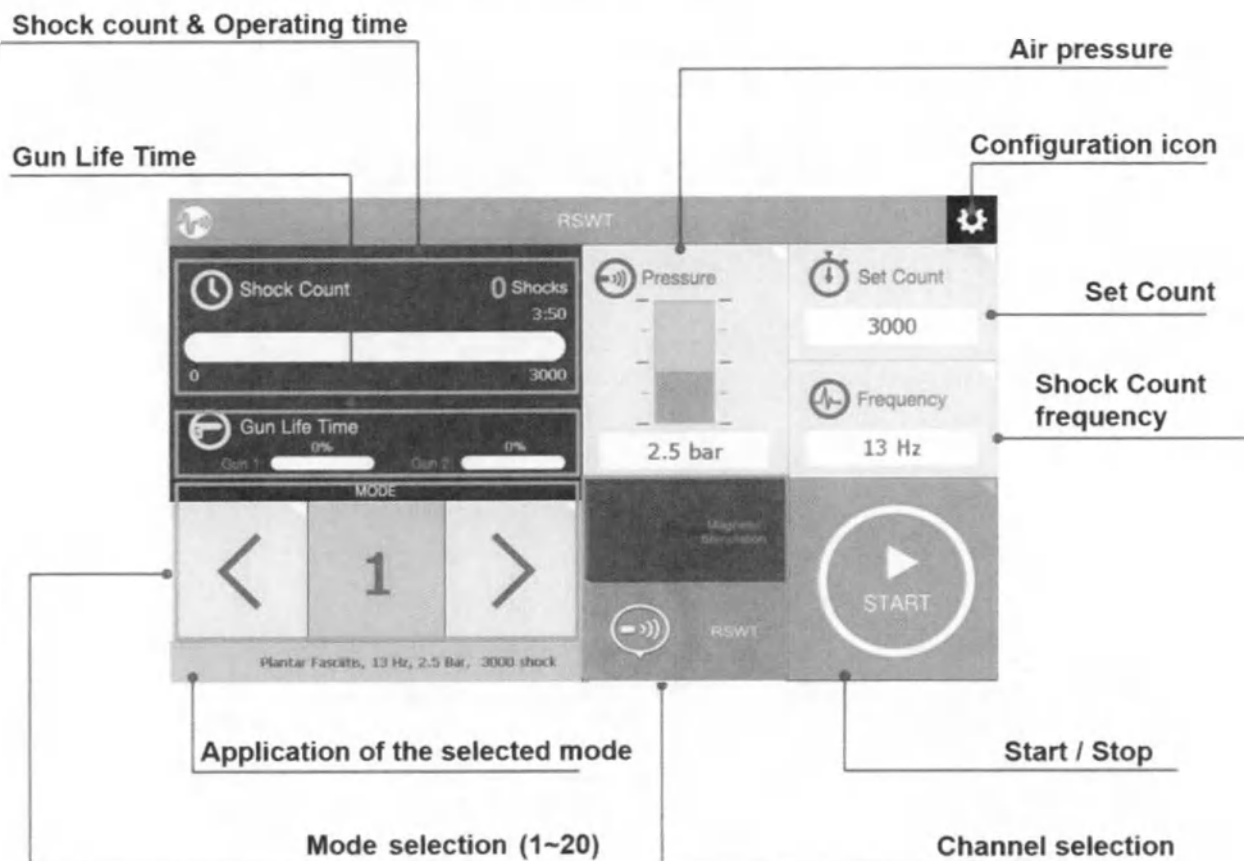
Figure 9. Controls - RSWT

5. Preparations for operation

5.1 Moving and fixing equipment

- ❑ Before moving the equipment, remove the connected power cable and the connection cables for peripheral equipment.
- ❑ When moving the equipment, place the transporting wheel at the lower section in the “loose” position.
- ❑ When the movement is completed, prevent the equipment from shaking by placing the transporting wheel in the “locked” position.
- ❑ If the equipment is moved while the wheel is in the “locked” position, it will damage the wheel, so please be careful.

5.2 Power connection

- ❑ Confirm that power socket and power cable are grounded before operating the machine.
- ❑ Plug the power cable into the power terminal of the OPTIMUS-Pro, and connect the other end to a grounded power socket.
- ❑ Confirm that the power cable is properly connected to the power terminal and power socket. If the connection is incorrect, unexpected problems can occur.
- ❑ Do not connect several devices to one power outlet.
- ❑ In particular, do not use a power cable that generates connection noise.

5. Operation instructions

6.1 Transducer placement

- You can adjust the position and angle of the transducer according to the patient and the stimulation area. The transducer arm supports can be moved up and down and left and right, and the angle of the transducer can be adjusted using the transducer holder. The height of the transducer arm is adjusted by the weight of the transducer using gas spring pressure. Refer to Figure 10 to adjust the height. To adjust the angle of the transducer, hold the transducer holder in the unlocked position and adjust the angle. The holder must be locked after adjustment is complete. Please refer to Figure 10 below.

Figure 10. Changing the angle of the transducer

- If the transducer hose is deflected as shown in Figure 11, the hose may break or be damaged. When positioning the transducer, make sure that the hose is not bent.

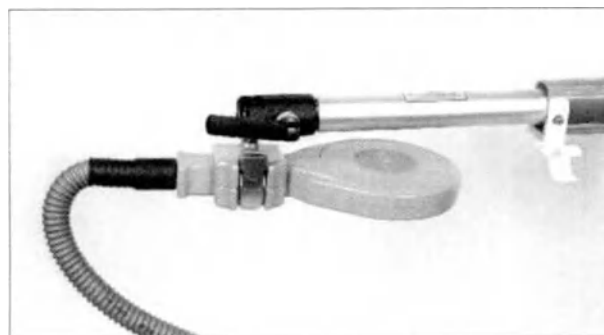
Figure 11. Deflected transducer

	Warning	Do not use the transducer if the transducer hose is bent.
--	----------------	---

6.2 Applications

The OPTIMUS-Pro can be used at different sites of the body depending on the pain. Figure 12-13 shows the use of the most common applications (Please note that there are some variations in the design of the transducer arm depending on the time of release.).

Figure 12. Application of the magnetic stimulation

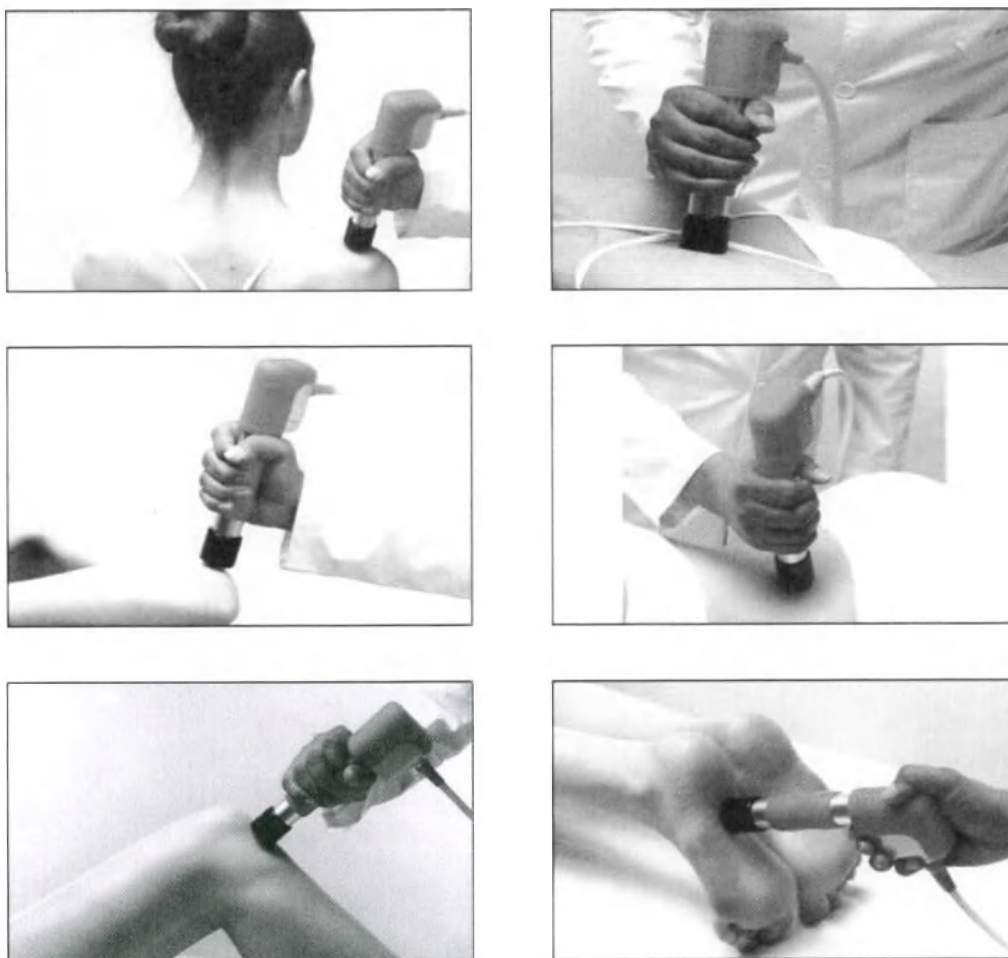
Figure 13. Application of the RSWT

6.3 Operation mode

Turn on the power switch located on the back of the equipment. When the power is turned on, the screen of the control panel is turned on as shown below, and the REMED logo screen appears for a moment, then it switches to the main screen.

Figure 14. REMED logo

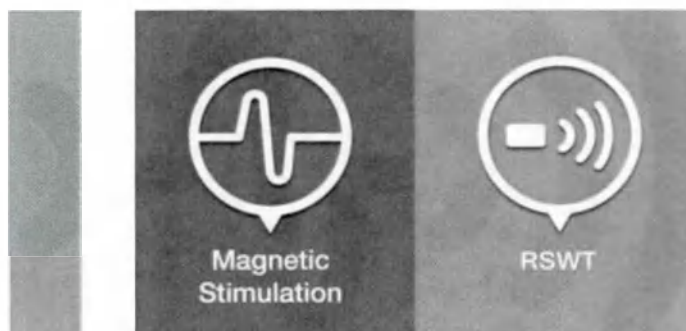
Figure 15. Function selection screen

Select the treatment function you want to use by pressing Magnetic Stimulation or RSWT on the function selection screen. Conversion of treatment functions is also possible on the main screen of each device operation, and both functions can be used at the same time

6.3.1 Magnetic stimulation function

Pressing Magnetic Stimulation on the initial function selection screen or RSWT function main screen switches the equipment to the electromagnetic function screen. The main screen configuration of the electromagnetic field function is as follows.

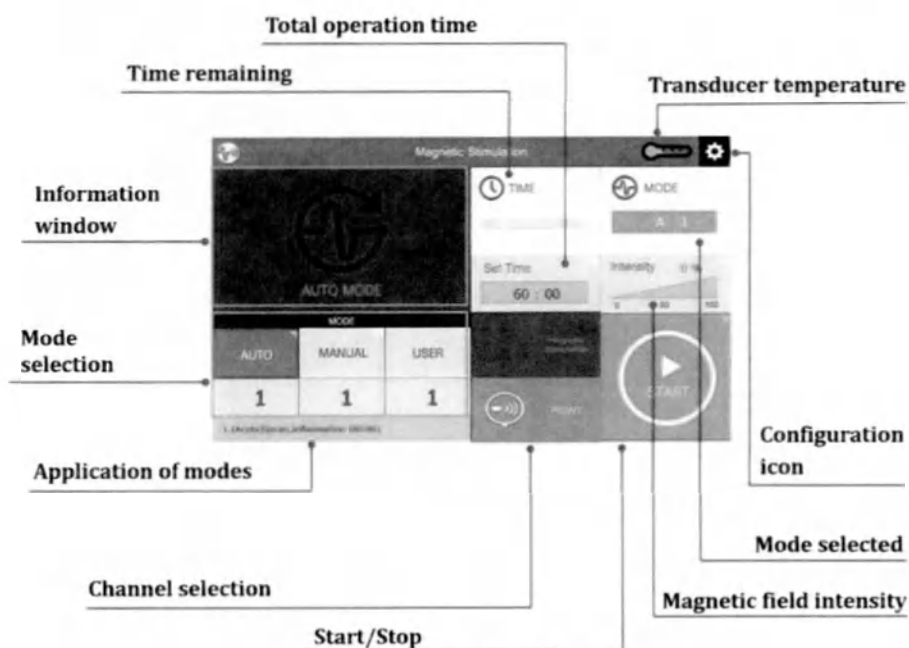
Figure 16. Magnetic stimulation main screen

The operating mode has Auto, Manual, User mode. While the Auto mode machine already sets the values of the parameter, users can set them in the Manual and User mode. Press Select to choose the mode.

(1) Auto mode

This mode has A1~A20 modes, and the frequency, train duration, inter train interval are set in advance. The waveform of selected mode is applied repeatedly for operating time by a user. Please see '8. Treatment protocols'. The ways of operation for Auto mode are like the following.

	(a) By using a touch screen, choose Auto mode and turn an encoder to choose the mode. The selected mode is presented upper and lower stage of a display..
	(b) Please set the time by turning the encoder to set the total operation timing. After press Start button, the equipment works at the intensity of magnetic field set by turning the encoder.
	(c) Intensity of the magnetic field can be also shown in the section of intensity.


(2) Manual mode

The Manual mode has M1 ~ M20 and the user is able to set directly the values of parameter. The parameter setting ranges are as follows.

- F1 : 1 ~ 100 Hz
- F2 : 1 ~ 100 Hz
- Ton : 1 ~ 4s (The setting ranges are different according to the value of F1 and F2.)
- Toff : 0 ~ 10s (The setting ranges are different according to the value of F1 and F2.)

The ways of instruction for Manual mode are like the following.

	(a) Please press the Manual mode for 2 ~ 3 seconds to set the parameters. on touch screen.
	(b) The user is able to set the parameter values after appearance of parameter setting screen.
	(c) Press Region of the mode and a screen will appear in which a name of the region can change. After changing the name, press Exit to close the screen.

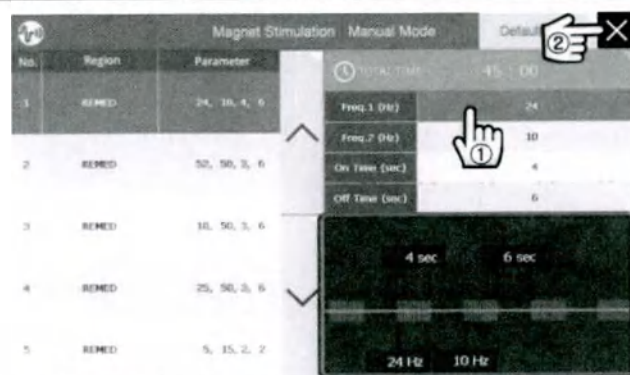(d) Choose the parameter to change. Freq1, Freq2, On Time, Off Time can be set in this mode. By turning an encoder, parameter values are set. Meaning of changed parameter value is shown in the bottom. After finishing all the settings, close the parameter setting screen by pressing a closing icon at the upper right corner.


(e) The operation method from this point is same as (b) ~ (d) in Auto mode.

(3) User mode


This has U1 ~ U20 mode, the user can set the parameter values directly. Each parameter setting ranges are as follows.

- Freq : 1 ~ 100 Hz
- On Time : 1 ~ 4s (The setting ranges are different according to the value of Freq.)
- Off Time : 0 ~ 10s (The setting ranges are different according to the value of Freq.)
- Strength : 1 ~ 100% (It is multiplying factor of the main magnetic field intensity.)
- Repetition : 1 ~ 100 cycle


The detailed installation and instruction of User mode are like the following.

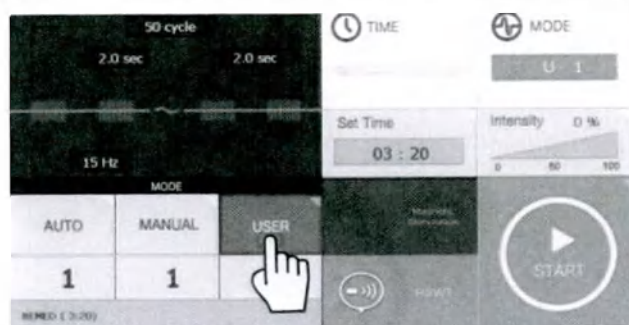
(a) In the beginning of using the equipment, press User mode on the touch screen for 2 ~ 3 seconds to set the parameter values.

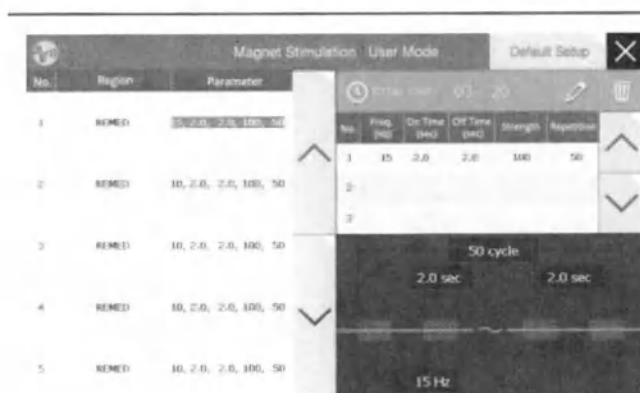
(b) The user is able to set the parameter values after showing parameter setting screen.

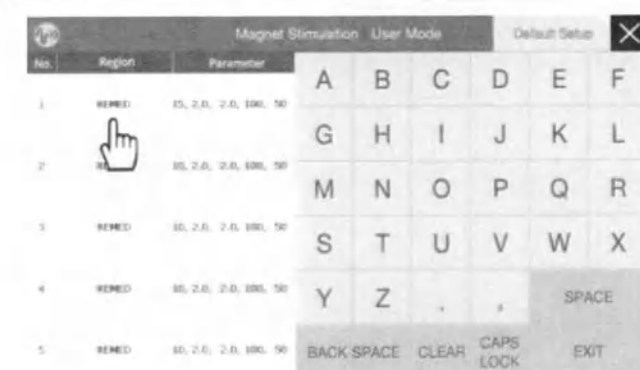
(c) By pressing Region of the mode for setting, alphabet screen appears and the name of the region can change.

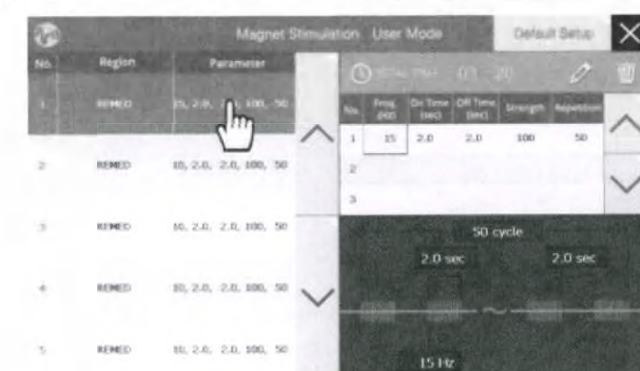
(d) By pressing the Parameter section, setting screen appears and the parameter values change.

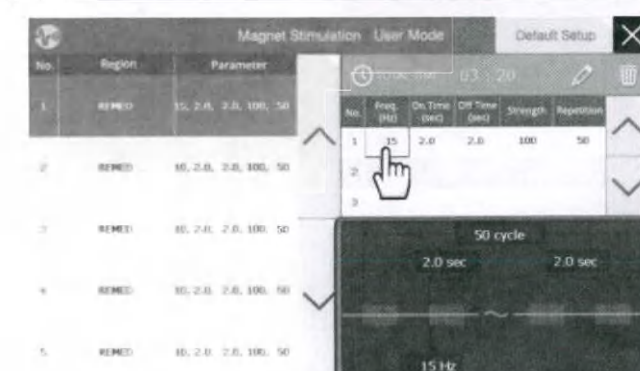(e) Choose the parameter to change. Freq1, On Time, Off Time, Strength and Repetition can be set in this mode. By turning an encoder, parameter values change. Meaning of changed parameter value are shown in the bottom.

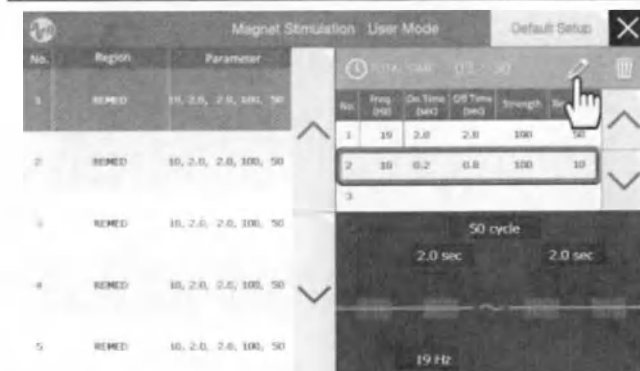(f) Each mode has 9 sub-modes with different parameter sets. The mode operates sequentially from sub mode 1. A sub-mode is added by pressing a pencil figure at the end of the top right parameter setting screen.

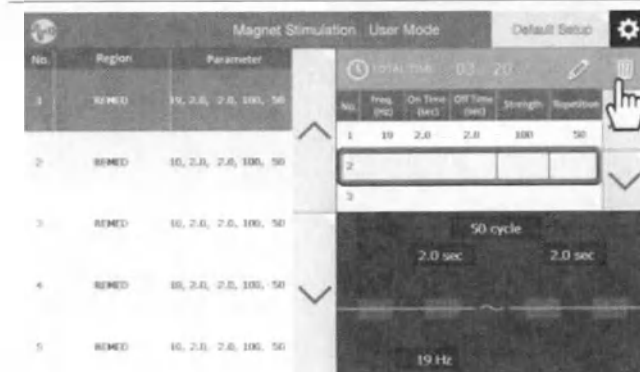(g) To remove the sub-mode, press a bin figure in the right top corner after choosing that sub-mode. Close the parameter setting screen by pressing a closing window in the end of the top right. After finishing all the settings, close the parameter setting screen by pressing a closing icon at the top end of the rightscreen.

(h) The operation method from this point is same as (b) ~ (d) in Auto mode.

6.3.2 Extracorporeal shockwave function

Pressing RSWT on the initial function selection screen or magnetic stimulation function main screen switches the equipment to the RSWT function screen. The main screen configuration of the RSWT function is as follows.

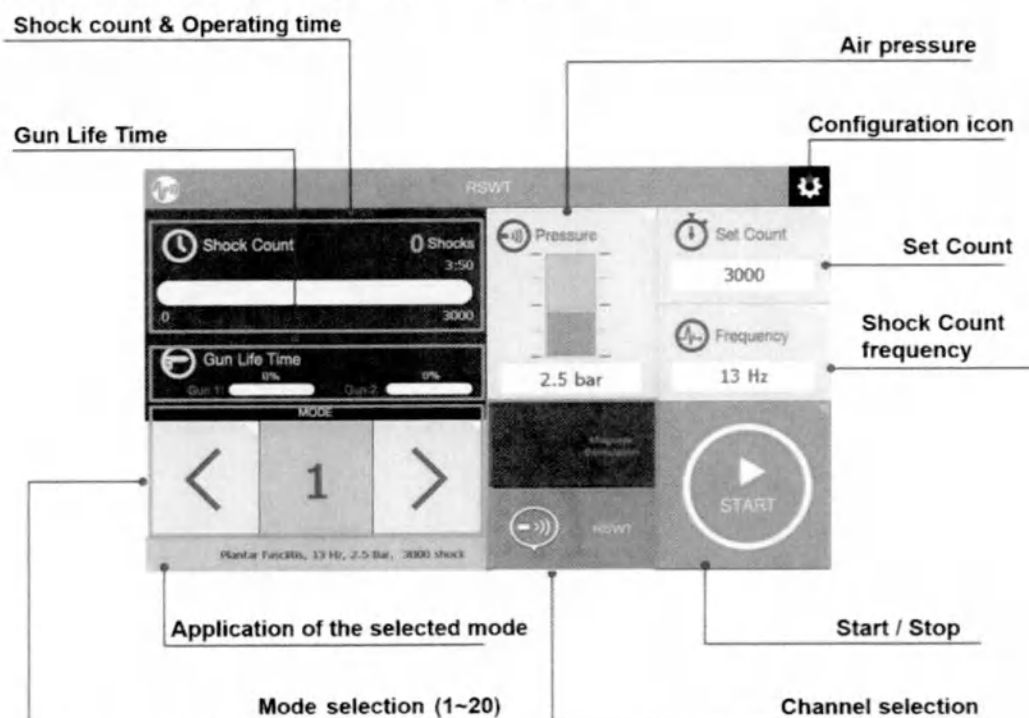
Figure 17. RSWT main screen

	Warning	Using the gun in the air without contact with the body may cause malfunction. To check the status of the gun, please use the lowest intensity of the shockwave when using in the air.
--	----------------	---

There are 2 channels and each has 20 modes. The user can set parameters, for example a number of shockwaves (Set Count), shockwave pressure (Pressure), Shockwave frequency (Frequency). The ranges of each parameter are as follows.

- Set Count : 100 ~ 5,000 shocks
- Pressure : 1 ~ 6 bar
- Frequency : 1 ~ 21 Hz

Operation instructions of each channel are same as follows and see '8. Treatment protocols'. The use of operation modes is like the following.

(a) Select the mode by using a left/right arrow or turning encoder. Application of the selected mode, parameter values and type of the applicator are shown.


(b) To temporarily change the parameter value, press the window of each variable to select it and turn the encoder to change the value.

However, when you turn the unit off and back on, or change to another channel, and then return to the original channel, the parameter value that was temporarily changed is not saved, so the value before the change appears.


(c) to (j) are methods to save parameters after changing parameters.


(c) When the user wants to change the parameter values, press mode number for 2 ~ 3 seconds.

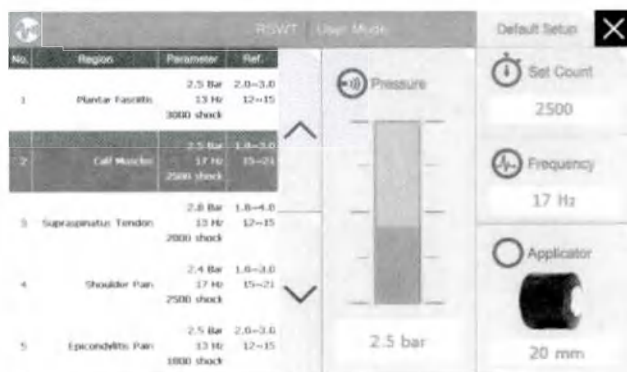
(d) The Parameter Settings screen appears, allowing the user to change parameter values.

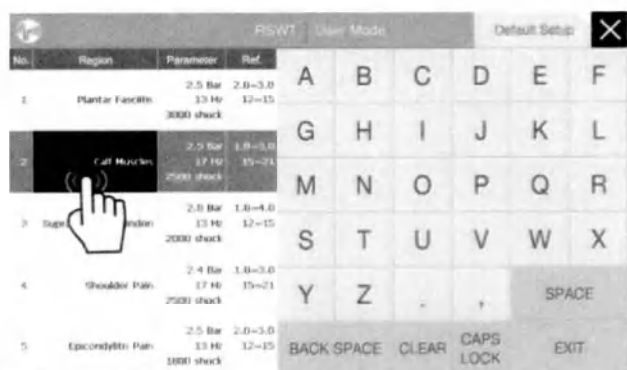
(e) Press Region of the mode for 2~3 seconds and a screen will appear in which a name of the region can change. After changing the name, press Exit to close the screen.

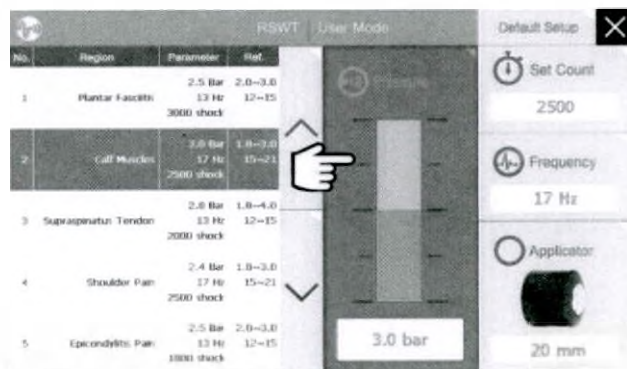(f) Set Count, Pressure, Frequency can change by selecting them and turning the encoder to change the values.

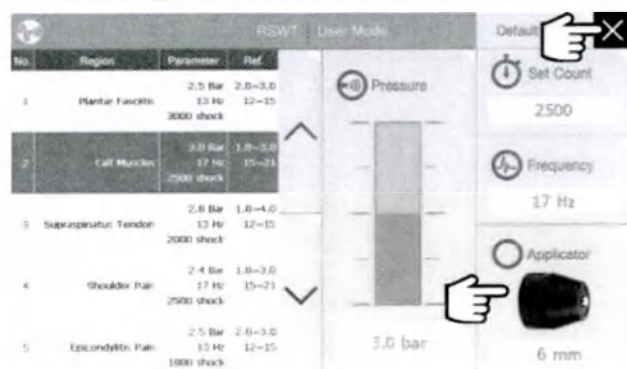(g) The user presses the applicator on the screen repeatedly until the type of the applicator required is displayed. When the user has finished setting the parameters, close the parameter setting screen by clicking closing icon at the upper right corner.

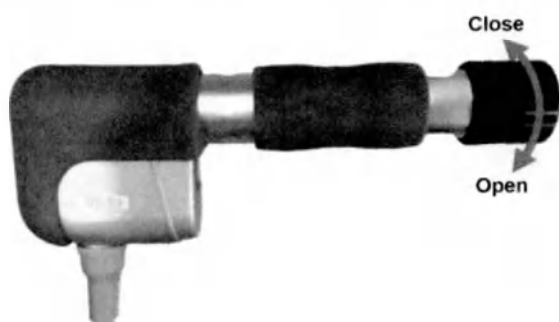(h) Disconnect the applicator connected with the Gun by turning clockwise the applicator. Install the applicator selected in (a) by turning anti-clockwise the new applicator.

(i) The changed parameter value and the selected applicator can be seen at the bottom of the mode window, in the parameter window, and in the applicator display of the channel selection.

(j) If you press the start button in the display once and then press the start button on the gun once, a shock wave will be generated (1 Shock). Press the start button on the gun for 2 to 3 seconds to start operation of the equipment. Pressure (Bar) and Frequency (Hz) can be changed during operation.

(l) To temporarily stop the operation of the equipment, press the start button on the gun once or press the encoder. If you want to restart the operation, press the start button on the gun again or press the encoder.

(m) Press the stop button on the touch screen when you want to stop completely before the equipment has finished its operation.

(n) When the operation is completed by the number of shockwaves generated, the operation of the equipment is terminated.

6.3.3 Environment setting

The environment setting of this equipment shows the number of gun output, the total number of gun output, the brightness of screen, and the intensity of input/output sound. You can reset the number of output generated from guns currently in use and adjust the brightness and sound. The setting method for this is as follows.

(a) Using the touch screen, press the configuration icon located at the top right of the screen.

(b) When the configuration screen appears, use the left and right indicators to adjust the brightness and sound of the screen. To reset the number of outputs generated from the gun currently in use, press the Reset button. When you are finished, click Close Window in the upper right of the screen to close the window.

7. Messages

7.1 Over Temperature - Transducer

The OPTIMUS-Pro is an equipment that generates a magnetic field by applying a high current to a transducer. As heat is generated from the transducer using the high current used to create the magnetic field, cooling is performed by circulating the cooling oil inside the transducer. If the transducer is used in a tilted position, the cooling oil will not reach all sections of the transducer, and injury or stoppage of operation can occur. Because of this, it is recommended to place the transducer in a horizontal position. When transducer is overheated, the operation is temporarily suspended and the "Over Temperature" message as shown in Figure 18 is displayed. When this message is displayed, put the transducer away from the patient for a while (Do not turn off the equipment) until the message disappears and the equipment returns to a normal state.

Figure 18. Over temperature messages – Transducer

	Caution	It is highly recommended to use the transducer in a horizontal position.
--	----------------	--

7.2 Transducer Error

The OPTIMUS-Pro always checks the connection of the transducer cable. If the cable is disconnected or damaged, the message is shown, Figure 19. In this case, equipment cannot be used anymore, and the user must contact authorized personnel of our local distributor maintenance.

Figure 19. Transducer error message

Note

The equipment always checks the connection state of the cable connected to the transducer.

7.3 Over Temperature - Air compressor (RSWT)

If the air compressor in the equipment is overheated during operation, a warning as shown in [Fig. 20] will appear and the generation of shockwave will be temporarily stopped. In this case, do not turn off the power of the equipment, wait for a while until the message disappears while the power is on, and it will be re-enabled. If you see this message frequently, please contact the REMED Customer Service Center for a checkup.

Figure 20. Over temperature - Air compressor (RSWT)

7.4 Gun Open

This equipment always checks the connection status of the gun. When the gun corresponding to the channel in use is not connected, the message "Gun Open" appears on the main screen as shown in Figure 21, 22 below.

Figure 21. Gun open message in RSWT function screen

Figure 22. Gun open message in magnetic stimulation function screen

7.5 Components exchange messages

Depending on the number of shocks that are output, this equipment shows you the message for replacement of components as shown in Figure 23, 24. Even if the message appears, the equipment works, but it is recommended to replace the components for safe use of the equipment.

(1) Tube & Projectile exchange

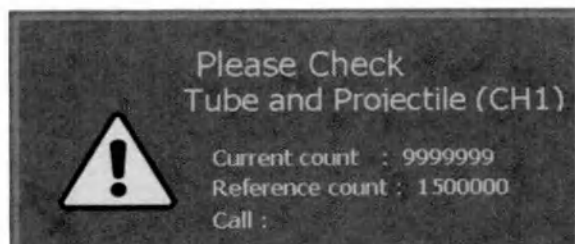
Figure 23. Tube & Projectile exchange message

(2) Gun exchange

Figure 24. Gun exchange message

J. Technical specifications

8.1 Size and weight

- ☐ Size (H x L x W): : 1960 x 515 x 580 mm
- ☐ Weight: : 76 kg

8.2 Power

- ☐ Input power : 230 VAC, 50 Hz
- ☐ Power consumption : Max 1.5 k VA

8.3 Modes – Magnetic stimulation

- ☐ Auto Mode : A1 ~ A20
- ☐ Manual Mode : M1 ~ M20
- ☐ User Mode : U1 ~ U20

8.4 Modes - RSWT

- ☐ Manual Mode : M1 ~ M20

8.5 RSWT number of shocks

- ☐ 100 ~ 5,000 shocks

8.6 RSWT intensity

- ☐ 1 ~ 6 bar

8.7 RSWT frequency

- ☐ 1 ~ 21 Hz

8.8 RSWT applicator

- ☐ 6, 10, 15, 20, 35 mm

8.9 Life Time

- 1 patient = 2000 shocks

1day = 3 patients = 6000 shocks

3year = 4,680,000 shocks

FSWT Probe Life Time = 5,000,000 shocks

= 3 Years

J. EMC emissions and immunity

Guidance and Manufacturer's Declaration – Electromagnetic Emissions

The OPTIMUS-Pro System is intended for use in the electromagnetic environment specified below. The customer or the user of the equipment should assure that it is used in such an environment.

The OPTIMUS-Pro requires special measures to ensure electromagnetic compatibility and must be installed and commissioned in accordance with the EMC-related information provided in this user manual.

Samples that require EMC compliance: power cord, shock wave gun with cable, display, holder arm with transducer holder assembly, electromagnetic transducer with hose.
! The use of transducers and cables that differ from those installed is not allowed. This can result in increased electromagnetic emissions or reduced noise immunity.
The OPTIMUS-Pro system should not be used in close proximity to or in conjunction with other equipment.

Emission Test	Compliance	Electromagnetic environment - guidance
RF emissions EN55011	Group 1	The OPTIMUS-Pro uses RF energy only for its internal function. Therefore, its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment.
RF emissions EN55011	Class A	The OPTIMUS-Pro system is suitable for use in all establishments other than domestic and those directly connected to the public low voltage power supply network that supplies buildings used for domestic purposes.
Harmonic Emissions IEC 61000-3-2	Class A	
Voltage fluctuations/ flicker emissions IEC 61000-3-3	Complies	

Guidance and Manufacturer's Declaration – Electromagnetic Immunity

The OPTIMUS-Pro System is intended for use in the electromagnetic environment specified below. The customer or the user of the equipment should assure that it is used in such an environment.

Immunity Test IEC 60601	IEC 60601 Test Level	Compliance Level	Electromagnetic environment - guidance
Electrostatic discharge (ESD) IEC 61000-4-2	± 8 kV contact discharge ± 15 kV air discharge	± 8 kV contact discharge ± 15 kV air discharge	Floors should be wood, concrete or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30%.
Electrical fast Transient/ burst IEC 61000-4-4	± 2 kV for mains cables ± 1 kV for input and output cables	± 2 kV for mains cables ± 1 kV for input and output cables	Mains power quality should be that of a typical commercial or hospital environment.
Surge	± 1 kV differential	± 1 kV differential	Anti-surge protection needs

IEC 61000-4-5	mode ± 2 kV common mode	mode ± 2 kV common mode	to be incorporated into the main supply to the equipment if surge protection is to be guaranteed.
Voltage dips, short interruptions and voltage variations on power supply input lines. IEC 61000-4-11	<5 % U_T (>95 % dip in U_T) for 0,5 cycle 40 % U_T (60 % dip in U_T) for 5 cycles 70 % U_T (30 % dip in U_T) for 25 cycles <5 % U_T (>95 % dip in U_T) for 5 sec	<5 % U_T (>95 % dip in U_T) for 0,5 cycle 40 % U_T (60 % dip in U_T) for 5 cycles 70 % U_T (30 % dip in U_T) for 25 cycles <5 % U_T (>95 % dip in U_T) for 5 sec	Mains power quality should be that of a typical commercial or hospital environment. If the user of the OPTIMUS-Pro System requires continued operation during power interruptions, it is recommended that the OPTIMUS-Pro System be powered from an interruptible power supply.
Power Frequency (50/60Hz) Magnetic Field IEC 61000-4-8	3 A/m	3 A/m	Power frequency magnetic fields should be at levels characteristic of a typical location in a typical commercial or hospital environment
Radiated RF Immunity EN61000-4-3	80 MHz – 2.5 GHz 2Hz 80% amplitude modulation		Equipment should only be used in the vicinity of other equipment compliant with EN60601-1-2.
Conducted RF Immunity EN61000-4-6	0.15 MHz – 80 MHz 2Hz 80% amplitude modulation		Equipment should only be used in the vicinity of other equipment compliant with EN60601-1-2.

NOTE U_T is the a.c. mains voltage prior to application of the test level

Guidance and Manufacturer's Declaration – Electromagnetic Immunity

The OPTIMUS-Pro System is intended for use in the electromagnetic environment specified below. The customer or the user of the equipment should assure that it is used in such an environment.

Immunity Test	IEC 60601 Test Level	Compliance Level	Electromagnetic environment - guidance
Conducted RF IEC 61000-4-6	3 V _{rms} 150 kHz – 80	3 V _{rms}	Portable and mobile RF communications equipment should be used no closer to any part of the OPTIMUS-Pro System, including cables, than the recommended separation distance calculated from the equation applicable to the frequency transmitter.
Radiated RF IEC 61000-4-3	3 V/m 80 MHz – 2.5	3 V/m	
			Recommended separation distance $d = \left[\frac{3.5}{\sqrt{P}} \right] V_1$ $d = \left[\frac{3.5}{\sqrt{P}} \right] E_1$ 80 MHz – 800 MHz $d = \left[\frac{7}{\sqrt{P}} \right] E_1$ 800 MHz – 2.5 GHz Where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer and d is the recommended separation distance in meters (m). Field Strengths from fixed RF transmitters, as determined by an electromagnetic site survey^a should be less than the compliance level in each frequency range^b.

			Interference may occur in the vicinity of equipment marked with the following symbol: 
--	--	--	--

NOTE 1 At 80 MHz and 800 MHz, the higher frequency range applies.

NOTE 2 These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

a. Field strengths from fixed transmitters, such as base stations for radio (cellular/ cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the location in which the OPTIMUS-Pro System is used exceeds the applicable RF compliance level above, the OPTIMUS-Pro System should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as re-orienting or relocating the OPTIMUS-Pro System

b. Over the frequency range 150 kHz to 80 MHz, field strengths should be less than [V1] V/m.

Guidance and Manufacturer's Declaration – Electromagnetic Immunity

The OPTIMUS-Pro System is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. The customer or the user of the OPTIMUS-Pro System can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters) and the OPTIMUS-Pro System as recommended below, according to the maximum output power of the communications

Rated maximum output power of transmitter, W	Separation distance according to frequency of transmitter		
	150 kHz - 80 MHz $d = \sqrt[3.5]{P} V_1$ Where $V_1 = 3$	80 MHz - 800 MHz $d = \sqrt[3.5]{P} E_1$ Where $E_1 = 3$	800 MHz - 2.5 GHz $d = \sqrt[7]{P} E_1$ Where $E_1 = 3$
0.01	0.117	0.117	0.233
0.1	0.369	0.369	0.738
1	1.167	1.167	2.333
10	3.689	3.689	7.379
100	11.667	11.667	23.333

For transmitters rated at a maximum output power not listed above, the recommended separation distance d in meters (m) can be estimated using the equation applicable to the frequency of the transmitter, where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer.

NOTE 1. At 80 MHz and 800 MHz, the separation distance for the higher frequency range applies.

NOTE 2. These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people

10. FAQ

Q1) Can it be used for the area of the body with artificial joints or plate pins?

A) Due to the nature of the magnetic field deep heat will be generated, which may cause deep-seated burns.

Q2) Is it the same effect to use on clothes?

A) One of the characteristics of a magnetic field is that it is permeable with a similar therapeutic effect. Therefore, it can be used on clothes.

Q3) Are there any precautions when using on clothes?

A) If there is a metallic decoration on the clothes, the magnetic field may cause burns on the surface of the skin. (e.g., bra strap rings, underwear metal ornaments, metal- containing fibers).

Q4) What part should be used with care?

A) If the transducer with the lever fixed is moved, the rubber ball inside will wear out and the transducer may come off. Therefore, release the lever when it is not in use

Q5) Is the applicator semi-permanent?

A) Applicator should be purchased as a consumable item. Please be careful because it can shorten the life time when you use it without contact with the affected part of the body.
B) Do not use it with very high intensity (pressure) because it may shorten the life time. Please use it with appropriate intensity.

Q6) When should the projectile be replaced?

A) It guarantees over 1,500,000 with good maintenance. It is recommended to clean it regularly.

Q7) Do I have to use a gel?

A) It is recommended to use it to deliver the shockwave without loss to the affected area of the body.

Q8) How do I maintain the equipment?

A) Projectile guide (pipe) cleaning is basically recommended to clean the oil at least once a week. If you use it a lot, we recommend cleaning the same day. Please refer to 13.3 Cleaning procedure.

1. Maintenance

11.1 Routine inspection of equipment

- ❑ The covering of power line of equipment, gun/transducer connecting line, etc. shall not be peeled off and internal lines shall not be exposed, and shall not be damaged by impact from outside.
- ❑ There shall be no trace of oil leakage from transducer.
- ❑ Wash the outside of equipment so that there is no foreign material.
- ❑ The button for equipment operation, etc. must not shake.
- ❑ The various parts attached to equipment must not shake.
- ❑ If any of the above occur, contact customer service for help

11.2 Safety inspection

- ❑ In order to ensure safe use, internal cleaning should be performed once per year by a person authorized by this company.
- ❑ In order to ensure safe use, be sure to check the equipment including internal components and output voltage from the person who has been given authorization from the company once per year.
- ❑ Please clean the transducer/gun before storing it
- ❑ When storing the product for a long period of time, be sure to check the product before using it.
- ❑ Please note the followings regarding storage conditions.
 - Keep out of water
 - Keep away from direct sunlight
 - Do not store near heaters
 - Avoid locations subject to excessive shock or vibration, exposure to chemicals or explosive gases.

11.3 Cleaning procedure

11.3.1 General cleaning

The OPTIMUS-Pro and its accessories can be cleaned in variable ways. Avoid any damage or contamination of the equipment using the methods recommended below. Before cleaning the unit, be sure to turn off the unit. Cleaning of the inside of the machine is only allowed to the person authorized by us.

- ❑ Periodically clean the equipment's exterior surface, gun and touch screen with a soft cloth moistened with alcohol. Do not use abrasives, lacquers, thinners, ethylene, or oxides, as these materials may cause permanent damage to the equipment.
- ❑ Do not submerge any section of equipment into liquid or detergent. In addition, no liquid must be allowed to enter the equipment or accessories.
- ❑ When the operation for one patient is completed, clean the applicator with alcohol and smooth cloth.
- ❑ If you want to use a disinfectant or a method not mentioned above, it may damage the device. Please note that in the event of damage to the equipment using unauthorized materials, no warranty of any kind is warranted.

11.3.2 Cleaning applicators

Apply the applicator after applying Transmission gel on skin due to the characteristics of this device. It is important to use a sufficient amount of Gel because air bubbles are formed and the pain is felt when shock waves are generated. For this reason, it is necessary to do regular cleaning once a week, since the gel may adhere to the applicator after prolonged use and may affect the shockwave generation

For cleaning, use running water to clean. If it is hardly fixed, it is effective to immerse the applicator in water for 2 to 3 hours and then wash it. After washing, remove the water using a soft dry cloth and dry it at room temperature for over 1 hour. Remove the applicator as shown in Figure 25 and clean the front metal part of the gun with a soft brush

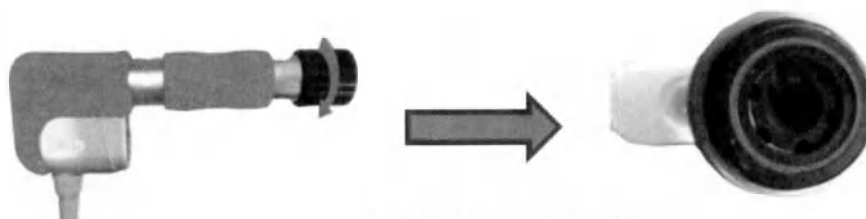
Figure 25. Applicator cleaning

(1) 6 mm Applicator cleaning

Disassemble the applicator from the gun as shown in Figure 26 below and clean the gel adhered to the arrow in the below.

Figure 26. 6mm Applicator cleaning

(2) 10 mm Applicator cleaning

Disassemble the applicator from the gun as shown in Figure 27 and clean the gel adhered to the arrow in the below.

Figure 27. 10mm Applicator cleaning

(3) 15 mm Applicator cleaning

Disassemble the applicator from the gun as shown in Figure 28 below and clean the gel adhered to the arrow in the below.

Figure 28. 15mm Applicator cleaning

(4) 20 mm Applicator cleaning

Disassemble the applicator from the gun as shown in Figure 29 and clean the gel adhered to the arrow in the below.

Figure 29. 20mm Applicator cleaning

(5) 35 mm Applicator cleaning

Disassemble the applicator from the gun as shown in Figure 30 below and clean the gel adhered to the arrow in the below.

Figure 30. 35mm Applicator cleaning

11.3.3 Projectile & projectile guide cleaning

As shown in Figure 31-(1) rotate the applicator to disassemble and rotate the center part of the gun to disassemble. Use the tweezers to remove the projectile guide in the center of the separated gun or the projectile left inside the gun, Figure 31-(2). Clean the projectile using a dry rag, and clean the inside of the gun using the cleaning tools provided Figure 31-(3).

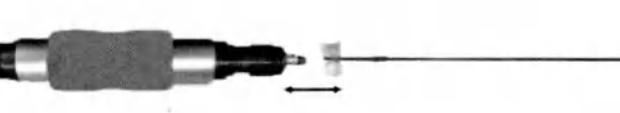
(1) Disassemble the applicator and gun.	
(2) Use the tweezers to remove the projectile guide in the center of the separated gun or the projectile left inside the gun. Clean the projectile using a dry rag.	
(3) Clean the inside of the gun using the cleaning tools provide.	

Figure 31. Projectile & projectile guide cleaning

11.3.4 Connection part of the gun cleaning

After a long period of use, a small amount of water and oil may flow out from the areas indicated by the arrows as shown in Figure 32. This is not a malfunction, so you can wipe it off with a soft cloth. However, if this happens frequently, please contact the REMED Customer Service Center for inspection.

Figure 32. Gun connection part cleaning

11.4 Replacement of components

This unit generates extracorporeal shockwaves using compressed air. Due to the nature of the machine, there are consumable components that need to be cleaned periodically and those that must be replaced due to long-term use. Details for exchange period and method are as follows.

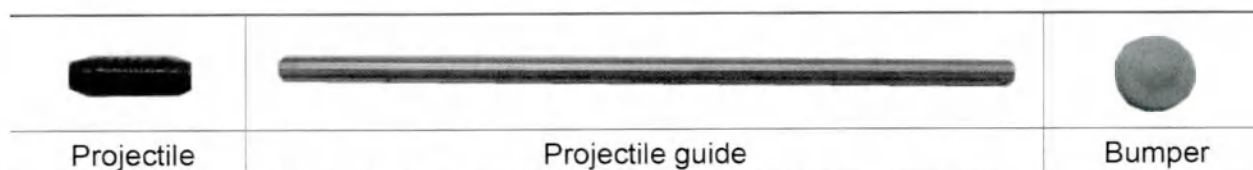
Figure 33. Projectile, projectile guide, bumper

(1) Exchange period

It is recommended to exchange after every 1,500,000 to 3,000,000 (3 bar, 15 Hz) use, but the exchange period may vary depending on the operation conditions.

(2) Exchange method

As shown in Figure 34-(1) rotate the applicator to disassemble and rotate the center of the gun to disassemble. Use the tweezers to remove the projectile guide in the center of the separated gun or the projectile left inside the gun, Figure 34-(2). The bumper can be also removed by using the tweezers. Use hex wrench to remove the projectile guide as shown in Figure 34-(3). After replacing new parts, assemble in reverse order.

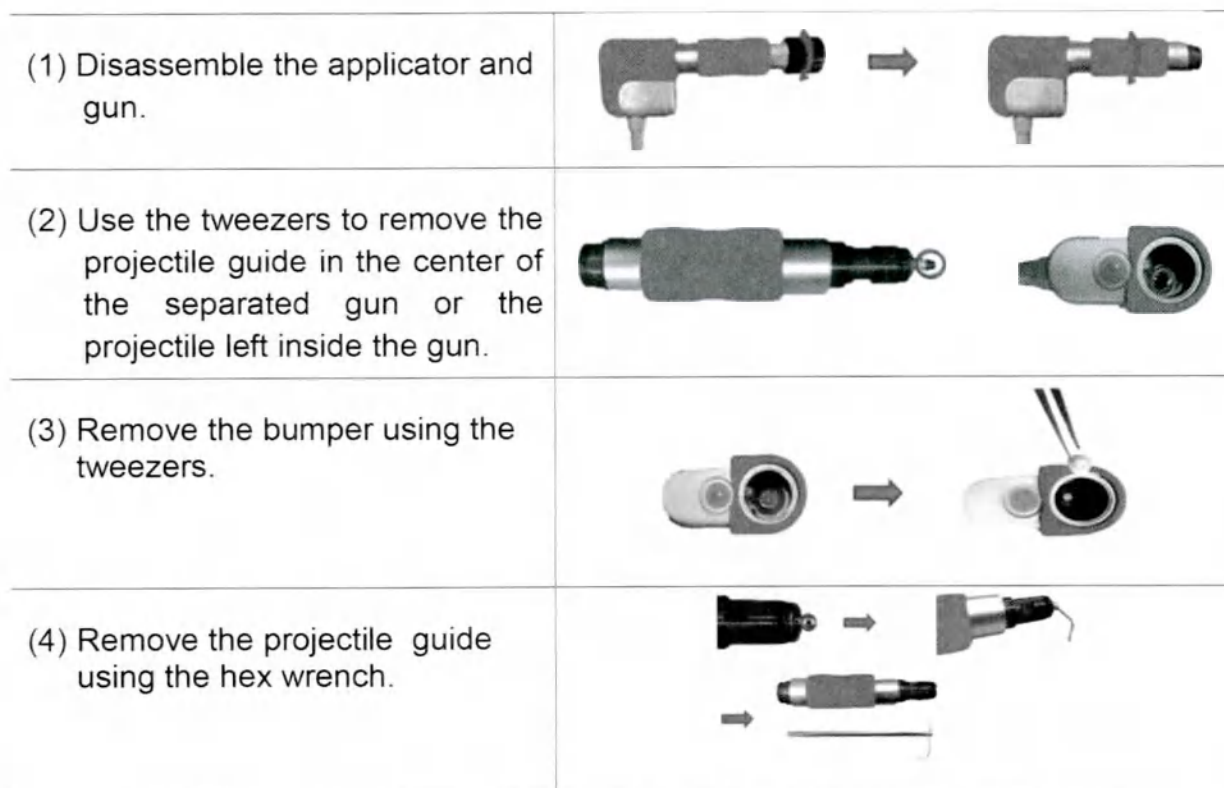
Figure 34. Projectile, projectile guide, bumper disassembling

11.5 Air compressor condensate remove

Air Compressor is an important component that compresses air and uses it in controlled amounts. The air stored in the air tank contains a large amount of moisture, oil, and dust. Therefore, the condensate in the air tank must be removed regularly

(3) Condensate remove period

Condensate should be removed after every 1,500,000 shocks of use. Otherwise, it may cause malfunction or malodor. Depending on the humidity of operating environment, condensate may need to be removed more frequently.

(4) Method of condensate remove

There is a black butterfly valve below the front of the equipment as shown in Figure 35. Turn on the power and turn the valve in the right (clockwise) direction while the pressure in the tank is raised

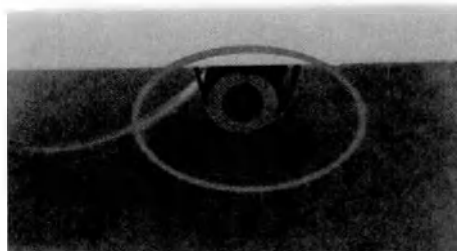
Figure 35. Air compressor condensate remove

When the condensate is removed, opening the valve can cause oil to come out and contaminate the environment. Combine the transparent hose provided with the valve and prepare the container to receive condensate. Then open the valve to remove condensate. When connecting the transparent hose to the valve, insert the hose into the blue one-touch fitting sleeve, Figure 36.

Please be careful because it may splash in all directions when condensate is removed. After all of the condensate is coming out, remove the hose by pressing the sleeve of the one-touch fitting as shown in Figure 37.

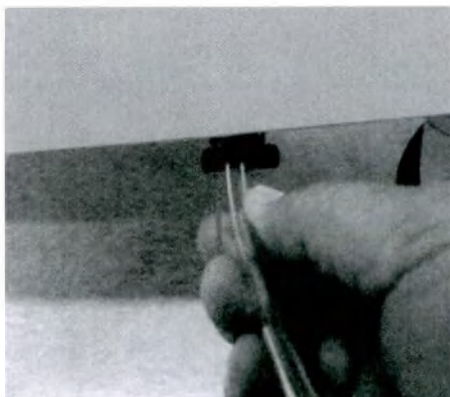
Figure 36. Hose connection

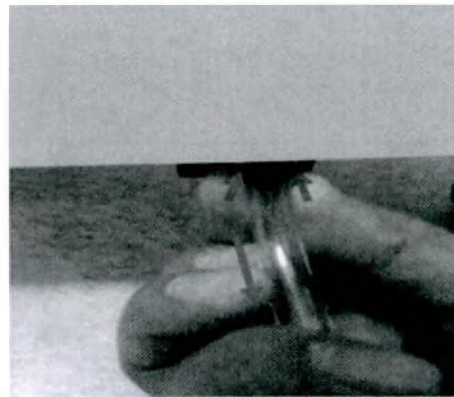
Figure 37. Hose disconnection

	Caution	After removing the condensate, remove the transparent hose from the one-touch fitting, and then turn the butterfly valve to the left (counter-clockwise) direction to lock it.
---	----------------	--

11.6 Troubleshooting

If the equipment does not operate normally during use, please check the items listed in the table below before requesting service. If none of the following problems apply, or if the following remedies do not help, turn off the power to the equipment and contact the REMED Customer Service Center.

Symptom	What to do	References in the Instructions For Use
Equipment does not turn on.	Check if the power connector of the equipment is properly connected.	• 3.4.2 Figure 1. Power cable connection
	Check if the power switch of equipment is turned on..	• 4.1.2 Figure 6. Rear view
	Confirm if the Key Lock switch ON.	• 4.1.1 Figure 5. Front view
	Confirm if the emergency button on the equipment is pressed. If it is pressed, turn it to the right to release it	• 4.1.1 Figure 5. Front view

Table 3. Self-troubleshooting

Symptom	What to do	References in the Instructions For Use
Magnetic field is not generated from equipment	Confirm if the LCD displays an ERROR message.	• 7. Messages
Over Temperature message is displayed.	Air compressor overheat message. Wait a moment and the message will disappear and be re-usable. If you receive this message frequently, please contact the REMED Customer Center for a checkup.	• 7.3 Over Temperature – Air compressor (RSWT)
	If transducer is in an upright position, it is easy for the transducer to overheat. Maintain a horizontal position as much as possible.	• 7. Messages
	If the room temperature is too high, disorder can be caused in the cooling. Maintain the room temperature at less than 25°C as much as possible.	• 3.3 Environment condition
When the gun is in operation, there is a windy sound.	Check the air tank valve lock. Turn the valve counterclockwise to lock it.	• 11.5 Figure 35. Air compressor condensate remove
	Please check the connection of gun connection cable..	• 4.1.2 Figure 6. Rear view
Gun's output is weak.	Moisture in the air tank can cause irritation. Please check and remove condensate.	• 11.5 Figure 35. Air compressor condensate remove
	If there is a dust inside the gun, or the applicator is loosened slightly, the output may be weak.	• 11.3.3 Projectile & projectile guide cleaning
There is no output from the gun.	Please check the connection of gun connection cable..	• 4.1.2 Figure 6. Rear view

Symptom	What to do	References in the Instructions For Use
	After you press Start on the touch screen, you must press the Start on the gun to get the output.	• 6.3.2 Extracorporeal shockwave function
	It may be relating the maintenance. (1) Replace the projectile guide and projectile. (2) Remove condensate (3) Check applicator assembly status Please.	• 11.4. Replacement of components • 11.5 Air compressor condensate remove
Wheel does not roll when moving the equipment.	Move the equipment by placing the wheel in the "Loose" position..	• 3.4.3 Figure 2. Transporting wheel in "Loose" and "Locked".

In the following cases, stop the operation by cutting off the power to the equipment, and contact the service center.

- The main power switch spontaneously turns off.
- The LCD screen of operation panel does not illuminate when power is turned off and then turned on again.
- When stimulation is not generated by gun.
- The temperature icon on the screen blinks and the equipment is not operated.
- The message indicating no connection between the gun and the main unit appears constantly.

11.7 Regular inspection for performance maintenance

In order to ensure safe use, be sure to do regular inspection from the person who has been given authorization from the company once per year.

11.8 Warranty

- ❑ This product is manufactured through an intensive quality control and inspection process.
- ❑ Compensation for repair and exchange of the product will be in accordance with the "Consumer Injury Compensation Rule" announced by the Economy Planning Board.
- ❑ The warranty period for this equipment is 1 year.
- ❑ If a failure occurs during the warranty period under normal operating conditions, the equipment will be repaired free of charge by our service center.
- ❑ If a problem with the equipment occurs during the warranty period, prepare the following customer service request to report to our company

Customer Service Request Form

The purpose of this form is to obtain information that will enable us to expeditiously process your service request. Please describe the fault or abnormality in as much detail as possible.

1. User information

- Name of hospital:
- Address of hospital:
- User name/telephone number: (Tel :)

2. Information on the equipment in use

- Name of model : OPTIMUS-Pro
- Date of purchase :
- Serial number of product :

3. Description of technical problems (Describe as much detail as possible.)

4. User checklist

If it is impossible to answer the question yes or no, please add information in the remarks column. If you are unsure of the answer, you can leave the item blank.

Number	Inspection item	Evaluation result	Remarks
1	Is use of the equipment impossible due to the current failure?	Yes/No	
2	When the power switch of the equipment is turned on, does the display function, or any part of the equipment function?	Yes/No	
3	Is output generated by the equipment?	Yes/No	
4	Is the strength of the output from equipment conspicuously different from when it was first installed? If so, how would you rank the difference in performance, in terms of a percentage?	Yes/No	
5	Can you perceive a change of output if the output strength of equipment is adjusted?	Yes/No	
6	Is there any leakage of liquid from equipment?	Yes/No	
7	Is there any damaged part in the external appearance? Describe this in the remarks column.	Yes/No	

After completely filling in the form, please return it to us.

Telephone : +(82) 42-934-5560

FAX : +(82) 42-934-5562

E-mail : service@remed.kr

12. Contact information

☐ Purchasing of product and technical inquiry

Address : 4th FL, SiliconPark A, 35, Pangyo-ro 255 beon-gil, Bundang-gu,
Seongnam-si, Gyeonggi-do, 13486, Republic of Korea

Telephone : +(82) 1588-7395

Fax : +(82) 31-696-4877

☐ Customer support

Address : #301~303, Migun Techno World II, 187, Techno 2-ro, Yuseong-gu,
Daejeon, 34025, Republic of Korea

Telephone : +(82) 42-934-5560

Fax : +(82) 42-934-5562

☐ European Representative

Obelis s.a.

Address : Boulevard Général Wahis 53, 1030 Brussels, BELGIUM

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This is hereby attested on this
26th day of Jul. 2023 at this office.

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Daehan Law Firm & Notary Office Inc.
Belong to Seoul Central District
Prosecutors' Office
DaeGwang B/D 3F, 36 Dongjak-daero.
Seocho-gu, Seoul, Korea

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
21st July 2023

Appendix to the User Manual

**Device for extracorporeal shock wave therapy (ESWT),
model SALUS-RSWT**

production: REMED Co., Ltd., Republic of Korea

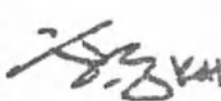Remed Co., Ltd.

Co-CEO
Eun-Hyun, Koh Sang-Yong, Lee

Appendix to the Instruction Manual
Extracorporeal shock wave therapy (ESWT) machine model SALUS
- RSWT

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WARNINGS AND PRECAUTIONS FOR USE:

- The installation and re-installation of the ESWT apparatus must always be carried out by specialized personnel authorized by the manufacturer.
- Attention! Make sure that the patient does not have any metal objects during the procedure. Failure to comply with this condition during the procedure can lead to serious skin burns!
- The installed ESWT apparatus is subject to regular safety checks by specialized personnel authorized by the manufacturer.
- Repair or installation of the ESWT apparatus should only be carried out by specialized personnel authorized by the manufacturer, and therefore any arbitrary disassembly/assembly of the equipment by the user is strictly prohibited.
- For electrical safety, the ESWT must always be connected to a properly grounded

power source in order to operate.

- The operator and manager of the ESWT apparatus must fully understand the instruction manual and it must be kept in close proximity to the equipment.
- To prevent safety-related accidents and to ensure proper maintenance, the safety signs, instructions, and regular inspection chart supplied with the ESWT apparatus should be posted in a conspicuous place near the ESWT apparatus.
- Since the ESWT has a fan inside the body to circulate air, remove curtains or any other items that may block the airflow around the ESWT.
- Do not use mobile phones, radios, mobile radio transmitters, cordless toys, etc. while using the ESWT device. in close proximity to the ESWT apparatus.
- During the operation of the ESWT device, the patient must not take drinks, water, etc., which may affect the operation of the ESWT device.
- Do not separate gun and applicator while ESWT is in operation, as this could result in injury due to air pressure release.
- If the transmission gel is not applied enough or if there is no contact with the skin, patients may feel pain because an air bubble is formed and a shock wave is generated. Therefore, it is necessary to apply a sufficient amount of gel and press the gun firmly against the patient to prevent pain.

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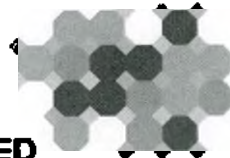
REMED

-
- in places where the temperature or humidity of the environment does not meet the requirements;
 - where corrosive or flammable gases may be present;
 - in places subject to strong shock or vibration;
 - in places directly exposed to wind and rain;
 - in dusty areas.

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Introduction	7
1 Purpose and contraindications	8
2 General description of the ESWT apparatus	11
3 Design features and appearance of the ESWT apparatus	12
4 Specifications	18
5 Materials used in manufacturing	20
6 Use of the ESWT device	21
6.1 Unpacking and installation	21
6.2 Preparing the ESWT apparatus for operation	22
6.2.1 Installation method	22
6.2.2 Power connection	22
6.2.3 Switching on	23
6.2.4 Shutdown	23
6.3 Operating mode parameters	24
7 Rules for care, cleaning and disinfection	25
8 Marking	26
9 Packaging	30
10 Safety requirements and environmental protection	31
1 1 Transport and storage	32
1 2 Maintenance and repair	33
1 3 Disposal information	35
1 4 Manufacturer's warranties	36
1 5 Claim information	38
1 6 Electromagnetic compatibility	39
List of regulatory documents used by the manufacturer to confirm compliance	41

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1 PURPOSE AND CONTRAINDICATIONS

Extracorporeal shock wave therapy (ESWT) machine model SALUS - RSWT **designed for** extracorporeal shock wave therapy . For more details, see "indications for use"

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Neurology:

- Pain syndromes: dorsalgia , cervicalgia , trigeminalgia , tension headaches, compression syndrome of the greater occipital nerve, piriformis muscle syndrome, intercostal neuralgia, tunnel syndromes.

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Lymphedema

Urology and gynecology:

- Overactive bladder
- Erectile dysfunction
- premature ejaculation

- Syndrome of acute and chronic pelvic pain
- Rehabilitation after prostate surgery
- Hypotension of the bladder
- Detrusor sphincter dyssynergy
- Hyposthesia of the glans penis
- Pudendal nerve pathology
- Stress urinary incontinence

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Decongestive therapy

Maxillofacial Surgery:

- Rehabilitation after CLS operations;
- shoulder syndrome;
- humeroscapular periarthritis ;
- physiotherapy for fractures ;
- chronic proximal plantar fasciitis in patients aged 18 years or older with symptoms for 6 months or more;
- plantar fasciitis ;
- traction degeneration of the plantar fascial band at the beginning of the medial tubercle of the calcaneus;
- repeated muscle contraction by passing electric currents through the electrodes without contact with the affected area of the body;
- removal of pain sensations (myofascial pain, thick fiber pain , fine fiber pain , neuropathic pain , phantom pain, fibromyogia);
- treatment of osteoporosis;
- treatment of fibromyalgia , arthritis, chronic fatigue.
- providing a completely non-invasive high-intensity pulsed electromagnetic effect on the musculature of the pelvic floor in order to rehabilitate weak pelvic muscles and restore neuromuscular control in the treatment of urinary incontinence in humans.

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- should not be used in conjunction with other electronic medical devices;
- **do not perform procedures in the area of the endoprosthesis;**

- The following patients must be approved by a doctor for the procedure:
 - patients with high body temperature (more than 37 °C), the elderly (over 85 years old) and children (under 3 years old);
 - person with implants (eg, artificial pacemaker, implantable infusion port, hearing aid).
- The following patients should not be treated with the product if:
 - blood disorders, clotting problems, or use of anticoagulants;
 - blood-thinning drugs (warfarin);
 - thrombosis;
 - oncological diseases;
 - with diabetes mellitus, diabetic neuropathy ;
 - pregnant women;
 - acute inflammation in the target area;
 - cartilage growth in children;
 - contact infectious patients (to avoid transmission of infection);
 - corticoid therapy for 6 weeks.

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2 GENERAL DESCRIPTION OF THE ESWT APPARATUS

The ESWT device consists of a body , power cord, shock wave therapy guns, and shock wave therapy applicators.

The body of the ESWT device is equipped with a key switch, an emergency stop button, an LCD touch screen (LCD -8") with a joystick , a compartment for storing applicators, rollers for moving with and without a brake, connectors for connecting shock wave therapy guns and a power cord, and as well as pre-installed software (SW).

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3 DESIGN FEATURES AND APPEARANCE OF THE ESWT DEVICE

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CH 1 (Gun 1) - gun shock wave therapy;

CH 2 (Gun 2) - gun shock wave therapy;

Emergency button – emergency stop button;

Controls – touch LCD display (LCD -8") with joystick;

Applicator box - a compartment for storing applicators;

CH 1 (Gun 1) connection port - connection connector of the shock wave therapy gun;

CH 2 (Gun 2) connection port - connection connector of the shock wave therapy gun;

power input terminal - power cord connector with holder bracket;

Figure 2. Power input terminal - power cord connector with holding bracket

power switch - power switch.

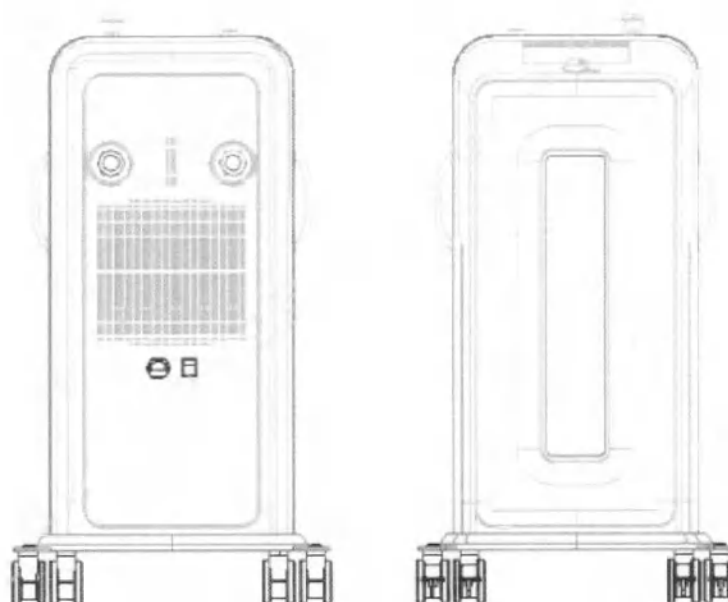
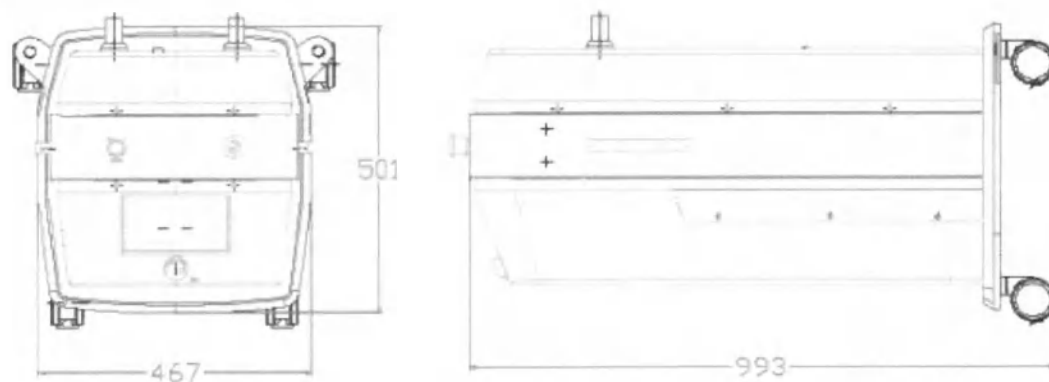
Figure 3. Dimensions of the body of the ESWT apparatus

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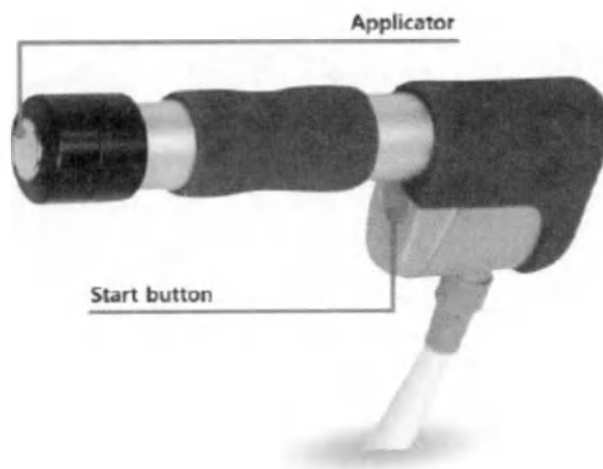
Figure 5. External view of the shock wave therapy gun

Applicator - applicator;
start button - start button.

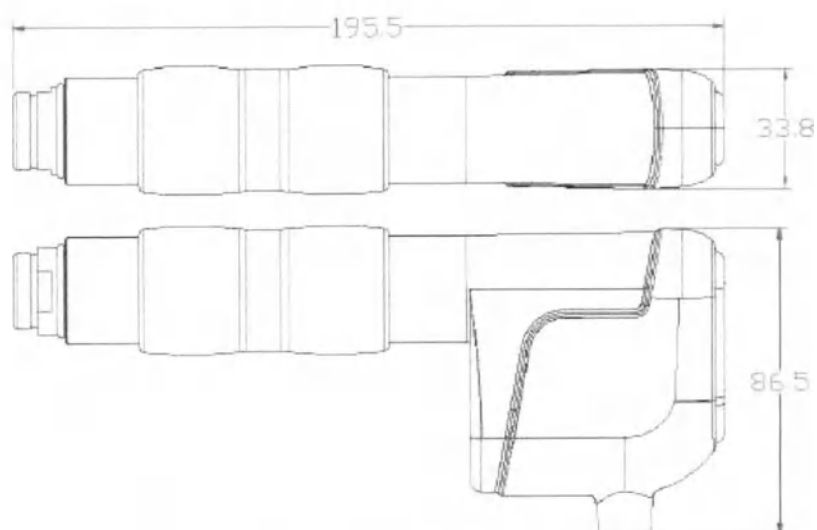
Figure 6. Shock Wave Gun Dimensions

6mm	10mm	15mm	20mm	35mm(option)
				

Figure 7. The appearance of the applicators and their size

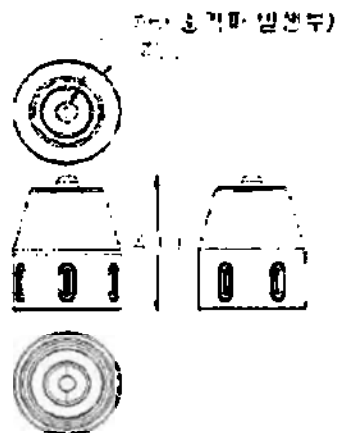
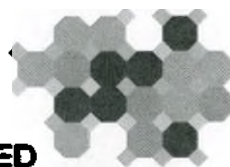
Figure 8. Applicator 6 mm

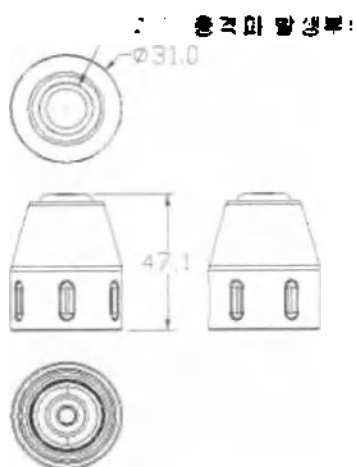
Figure 9. Applicator 10 mm

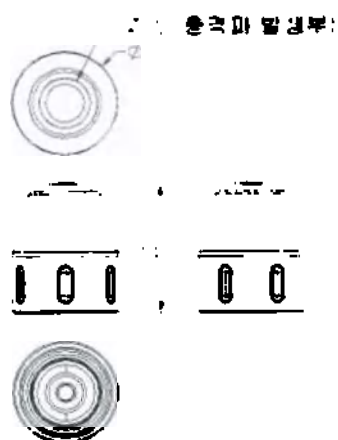
Figure 10. Applicator 15 mm

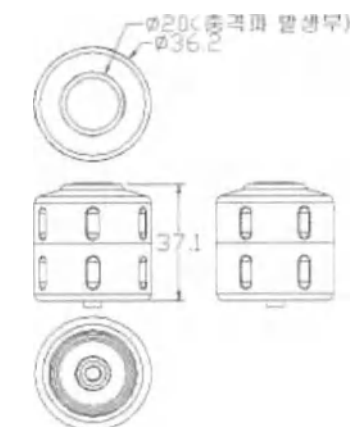
Figure 11 Applicator 20 mm

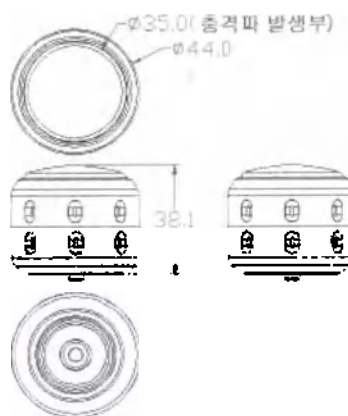
Figure 12 Applicator 35 mm

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Table 1: Software details

Name of the software	Software identification name	Version number (identification number) of the software	Digital software identifier (checksum of the executable code)	Software Digital ID Calculation Algorithm	Software release date
RSWT F/W (firmware firmware)	-	1.00.11	-	-	06/30/2017
RSWT S/W (software)	-	2.00.09	-	-	09/18/2017

Note: Firmware firmware and software perform different functions together in one product. Firmware is designed to control the equipment - controls the output according to the input from the software), the software is designed to control only what is shown on the display: therapy time, buttons.

 Warning

- **Software updates should only be carried out by an authorized representative. The software is encrypted and cannot be modified by yourself.**
-

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Table 2

Parameter name	Value, norm
1	2
Overall dimensions of the ESWT apparatus (H×L×W), mm	993 x 467 x 501
ESWT body dimensions (H×L×W), mm	993×467×501
Dimensions of the compartment for storing applicators (H×L×W), mm	70×242×72
Roller diameter with brake, mm	75
Number of rollers with brake, pcs	2
Roller diameter without brake, mm	75
Number of rollers without brake, pcs	2
Mass of the body of the ESWT apparatus, kg	56
Shock wave therapy gun dimensions (H×L×W), mm	33.8×195.5×86.5
Shock wave therapy gun cable length, mm	2000
Shock wave therapy gun holder dimensions (H×L×W), mm	49 × 77 × 10 6.4 _
Number of shock wave therapy guns, pcs	2
Weight of shock wave therapy gun, kg	0.6
Number of operating modes, pcs	1 (shock wave therapy)
Time of readiness for work after inclusion, min., no more	10
Maximum time of continuous operation, min., not less than	Long-term with SWT
Power cord length, mm	1500
Power cord weight, g	130
Rated supply voltage, V	220
Power consumption, kVA , no more	1.5
Nominal value of the mains frequency, Hz	50
Display type	built-in
Product weight without packaging / with packaging, kg	50/68
Joystick dimensions (Diameter×Height), mm	34×14
Emergency stop button dimensions (Diameter×Height), mm	40×28
Key switch base dimensions (Diameter×Height), mm	30×11
Key switch dimensions (Diameter×Height), mm	21×14
<i>Display specifications:</i>	
Display type	color TFT LCD (Touch Screen)
Dimensions (L×W), mm	180×103

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Table 3 Characteristics of applicators

Applicator size	Dimensions (height×diameter), mm	Impact surface diameter, mm	Weight, g
1	2		3
Applicator 6 mm	40.0×31.0	6	31
Applicator 10 mm	47.1×31.0	10	42.5
Applicator 15 mm	38.7×31.0	15	57
Applicator 20 mm	37.1×36.2	20	68
Applicator 35 mm	38.1×44.0	35	100.5

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Protection class against electric shock: I, working part Type B F.

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Service life - 3 years.

5 MATERIALS USED IN MANUFACTURING

The ESWT apparatus is made of the following materials:

- ABS plastic PA 200 – for body and body elements, shock wave therapy gun, applicators;
- Stainless steel 420 SS - for the body and body elements;
- XINJIANG grade SG-3 polyvinyl chloride for power cord, gun cable and transducer hose;
- Rubber grade IR -307 - for inserting the shock wave therapy gun and applicators;
- Stainless steel grade AMS 5355 by Wire Company , Inc. _ - for applicators.

The ESWT device does not contain medicinal products, components of animal or human origin.

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6.2.1 Installation method

- When moving the ESWT apparatus, rollers on the body are used;
- To move, set the braked castors to the open position (see Figure 4);
- After moving, set the brake rollers to the closed position (see Figure 4);
- Place the ESWT machine on a flat and level surface;
- Install the ESWT in a place with suitable ambient temperature and humidity, and do not install in a place where it will be exposed to dust or flammable materials. Operate at temperatures from +10°C to +30°C and relative air humidity from 30 to 85%;
- Connect the power cord to the body of the ESWT (see Figure 2);
- Connect the power cord to a grounded socket;
- Make sure the power cord is connected;
- Before turning on the equipment, make sure that the socket and power cable are grounded;
- Be careful not to damage the equipment with strong impacts;
- Plug in the shock wave therapy gun;
- Turn on the power switch located on the back of the ESWT (see Figure 1);
- When the power is turned on, the manufacturer's logo will appear on the display of the ESWT device for a moment and then the "Main Screen" appears;
- The patient is seated or placed in front of the apparatus;
- Connect the appropriate applicator to the shock wave therapy gun;
- The applicator is placed over the desired area of the body depending on the pain;
- Select the necessary parameters for the procedure;
- During the procedure, the display shows data on the parameters of shock wave therapy, as well as warnings;
- The duration of the procedure depends on the selected parameters.

Power connection

Please use the original power cord.

Figure 14. Power cord

Please make sure the AC power supply meets the following specification: AC 220(V), 50(Hz).

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The parameters of the operating modes correspond to table 4.

Table 4

Parameter name	Value, norm
1	2
Operating modes	Manual with the possibility of entering into the memory of the device 20 programs with adjusted parameters of shock wave therapy (from M 1 to M 20). Adjustable parameters: number of strokes, frequency, pressure
Frequency range, Hz with resolution, Hz and tolerance, %	1-21 1 ± 20
Pressure range, bar with resolution, bar and tolerance, %	1-6 0 ± 20
Range of number of beats, beats with discreteness, strokes	100-5000 100
Standard sizes of used applicators	6, 10, 15, 20, 35

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Warning

- Before cleaning the guns, turn off the ESWT machine and disconnect the gun cables from the connectors;
 - Preventive checks and cleaning of applicators, guns and external surfaces of the ESWT apparatus can be carried out by users;
 - Cleaning of the inside of the ESWT body is only permitted by personnel authorized by the manufacturer.
-

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Disinfectant - 70% isopropyl alcohol.

Warning

- NEVER immerse the ESWT device and its parts in water or any solution, pour water or any solution on the ESWT device and its parts;
 - NEVER allow liquid to enter the housing, switches, connectors, or any vents on the ESWT;
 - ALWAYS wipe off excess liquid remaining on the surface of the ESWT apparatus and its parts with a cloth;
 - NEVER use abrasive, corrosive cleaners or cleaners containing acetone;
 - NEVER disinfect the ESWT device and its components using high pressure and high temperature.
-

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
REMED

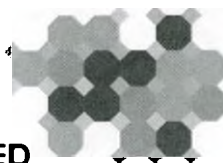8 MARKING**8.1 Marking applied to the label affixed to the body of the ESWT apparatus:**

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

8.2 Marking applied to the packaging of the device:

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
8.3 Decoding of all marking symbols:

Symbol	Meaning	Description
	Warning	Warning of danger to the life of the operator if the precautions are not observed
	Attention	Notice of precautions to be taken to avoid injury or damage
	Don't push	It is forbidden to push
	Don't sit	It is forbidden to sit on the surface
	Type BF	Working part type BF
	Work instructions	You must follow the operating instructions
	Operating instructions	You must follow the operating instructions
	Serial number	Product serial number

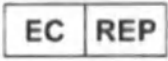Symbol	Meaning	Description
	Authorized Representative in the European Community	Authorized Representative in the European Community
	vendor code	Product number
	Manufacturer	Manufacturer
	date of manufacture	date of manufacture
	Do not stack	Note on the prohibition of stacking products in packaging
	Keep away from moisture	Instruction on the need to store the transport packaging away from rain and in dry conditions
	fragile	Fragile. Carefully
	Top	Specifying the Correct Vertical Position of the Shipping Package
	Reuse	Indication of the need for disposal at the appropriate collection points for waste to be recovered or recycled
	Off (power)	Designation of the position of the mains switch corresponding to disconnection from the mains
	Included (food)	Designation of the position of the mains switch corresponding to the inclusion in the network
	Alternating current	Indication on the product label about the need to connect the product through the appropriate connector only to AC sources
	Do not dispose of with household waste	Indication of the need for disposal at separate collection points for recovery and recycling

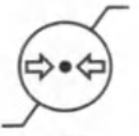Symbol	Meaning	Description
	Temperature Range	Indication of the maximum and minimum temperatures for storage, transport and operation
	Humidity range	Specify the permissible upper and lower limits of relative humidity during transportation and storage
	Atmospheric pressure range	Indication of permissible upper and lower limits of atmospheric pressure during transportation and storage
	dangerous voltage	Caution, dangerous voltage
	CE marking	CE marking

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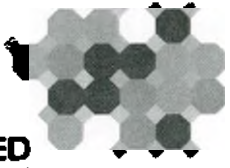The safety of the ESWT device complies with the requirements of IEC 60601-1.

For electromagnetic compatibility, the ESWT device complies with the requirements of IEC 60601-1-2 (class A, group 1).

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⚠ Warning

- The user should routinely inspect the ESWT device, cables, guns, and applicators for damage that could affect patient safety or performance. The recommended inspection interval is one week or less. Repair damage before using the instrument;
- Periodic safety checks must be carried out by personnel authorized by the manufacturer. The recommended interval is one year, or you must follow the hospital's established intervals and procedures.

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- The outer sheath of the ESWT power cord, gun cables, etc. must not peel off, the internal wires must not be exposed and must not be damaged by external impacts;
- Clean the outside of the ESWT apparatus so that there are no foreign objects;
- The joystick for controlling the ESWT apparatus, etc., must not tremble;
- The various parts attached to the ESWT apparatus must not shake.

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
35

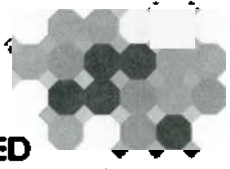REMED

In case of post-warranty and out-of-warranty repair of the ESWT device, contact the authorized representative of the manufacturer in the territory of the Russian Federation: Infomed-Neuro Limited Liability Company (Infomed-Neuro LLC) at the address: Russia, 125445, Moscow, Leningradskoe shosse, bld. 69, building 1, room I, fl . 1 room 50; phone +7 495 645 47 00; email kitayeva @ im - neuro . ru .

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Infomed-Neiro Limited Liability Company (Infomed-Neiro LLC) at the address: 1 Russia, 125445, Moscow, Leningradskoe shosse, 69, building 1, pom. I, fl . 1 room 50; phone +7 495 645 47 00; email kitayeva @ im - neuro . ru .

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REMED

15 COMPLAINT DETAILS

In case of warranty and out-of-warranty repair of the product, contact: Limited Liability Company "Infomed-Neuro" (LLC "Infomed-Neuro") at the address: Russia, 125445, Moscow, Leningradskoe shosse, 69, building 1, pom. 1, fl. 1 room 50; phone +7 495 645 47 00; email kitayeva @ im - neuro . ru .

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Warning

ESWT devices require special measures to ensure **ELECTROMAGNETIC COMPATIBILITY** and must be installed and put into service in accordance with the EMC-related information given in the instruction manual.

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
FAX : (02) 3477-0957

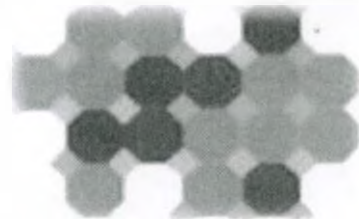REMEDI

Rehabilitation Equipment Company

APPROVE

CEO

REMEDI CO., Ltd.

Eun Hyun Koh, Sang Yong Lee

21st July 2023

Appendix to the User Manual

**Device for extracorporeal shock wave therapy (ESWT),
model SALUS - Optimus – Pro**

production: REMEDI Co., Ltd., Republic of Korea

Remed Co., Ltd.

Co-CEO
Eun-Hyun, Koh

Sang-Yong, Lee

Appendix to the Instruction Manual
Extracorporeal shock wave therapy (ESWT) machine model SALUS
- Optimus - Pro

(hereinafter referred to as ESWT device(s), product(s), device(s))

Before using this appliance, please read the instruction manual and this appendix to the manual carefully and save it for future reference!

SAFETY INSTRUCTIONS

ELECTRICAL SAFETY PRECAUTIONS:

- The rated voltage of the power supply for the ESWT apparatus is 220 V, alternating current;
- Make sure that all connecting parts (power cord or component cables) are properly connected to the ESWT machine;
- Make sure the ESWT apparatus is fully grounded;
- Repair, expansion and installation of the ESWT apparatus must not be carried out by anyone other than specialized personnel authorized by the manufacturer. Arbitrary disassembly / assembly of the ESWT device by the user is strictly prohibited;
- Before connecting other equipment not listed in the instruction manual, be sure to consult the manufacturer or an authorized representative of the manufacturer;
- To avoid electrical interference during use, the ESWT device should be installed at a considerable distance from the power source, any generator, x-ray equipment, transmitter, mobile electrical wire, etc.;
- Be careful of possible electromagnetic or other interference between the ESWT machine and other devices;
- As such, an independent power circuit is required and power sharing with other electronic devices is not recommended;
- Do not place the ESWT in a location where it is difficult to disconnect the power plug or other detachable plug.

READ THE OPERATING MANUAL AND THIS APPLICATION BEFORE USING THE APPLIANCE

WARNINGS AND PRECAUTIONS FOR USE:

- The installation and re-installation of the ESWT apparatus must always be carried out by specialized personnel authorized by the manufacturer.
- Attention! Make sure that the patient does not have any metal objects during the procedure. Failure to comply with this condition during the procedure can lead to serious skin burns!
- The installed ESWT apparatus is subject to regular safety checks by specialized personnel authorized by the manufacturer.
- Repair or installation of the ESWT apparatus should only be carried out by specialized personnel authorized by the manufacturer, and therefore any arbitrary disassembly/assembly of the equipment by the user is strictly prohibited.
- For electrical safety, the ESWT must always be connected to a properly grounded power source in order to operate.
- The operator and manager of the ESWT apparatus must fully understand the

instruction manual and it must be kept in close proximity to the equipment.

- To prevent safety-related accidents and to ensure proper maintenance, the safety signs, instructions, and regular inspection chart supplied with the ESWT apparatus should be posted in a conspicuous place near the ESWT apparatus.
- Since the ESWT has a fan inside the body to circulate air, remove curtains or any other items that may block the airflow around the ESWT.
- Since a strong magnetic field is generated around the magnetic field generation section, technicians, assistants, and patients should not hold any objects that may be affected by the magnetic field in their hands.

Note: Magnetic fields can damage items such as watches and mobile phones, so be careful and store them separately.

- Do not use mobile phones, radios, mobile radio transmitters, cordless toys, etc. while using the ESWT device. in close proximity to the ESWT apparatus.
- During the operation of the ESWT device, the patient must not take drinks, water, etc., which may affect the operation of the ESWT device.
- Do not separate gun and applicator while ESWT is in operation, as this could result in injury due to air pressure release.
- If the transmission gel is not applied enough or if there is no contact with the skin, patients may feel pain because an air bubble is formed and a shock wave is generated. Therefore, it is necessary to apply a sufficient amount of gel and press the gun firmly against the patient to prevent pain.

(Gels registered in the territory of the Russian Federation* can be used as a transmission gel (hereinafter referred to as the gel):

1) Ultrasound gel "Mediagel" with Aloe Vera juice according to TU 21.20.23-001-76063983-2019, registration certificate No. RZN 2021/15126 dated 19.08.2021.

2) Medical conductive abrasive gel "Akugel-Preparation" according to TU 32.50.50-013-66242897-2020, registration certificate RZN 2021/14167 dated 04/28/2021.

*** Not included in delivery.**


- The shock wave can be subjected to excess pressure when used around the lungs, ears, and gastrointestinal tract.
- Do not apply shock waves directly to the heart, head, or injured area.
- The ESWT machine should be used for patients who can accurately express and convey pain during the procedure.
- Excessive pressure on the treated area with a gun can cause skin damage such as blisters and redness.
- Do not move the ESWT device during use because the gun may fall or the ESWT device may be damaged due to impact or movement of the ESWT device.

Do not use the ESWT device in any of the following locations:

- in places where the temperature or humidity of the environment does not meet the requirements;
- where corrosive or flammable gases may be present;
- in places subject to strong shock or vibration;
- in places directly exposed to wind and rain;
- in dusty areas.

When using the ESWT in the following locations, use appropriate and sufficient protective measures:

- in places subject to static electricity or any other interference;
- in places exposed to intense electric or magnetic fields;
- in places located near power lines.

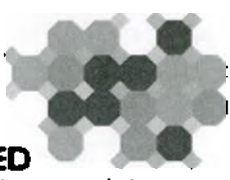REMED

CONTENT

Introduction	7
1 Purpose and contraindications	8
2 General description of the ESWT apparatus	12
3 Design features and appearance of the ESWT apparatus	13
4 Specifications	20
5 Materials used in manufacturing	23
6 Use of the ESWT device	24
6.1 Unpacking and installation	24
6.2 Preparing the ESWT apparatus for operation	25
6.2.1 <i>Installation method</i>	25
6.2.2 <i>Power connection</i>	26
6.2.3 <i>Switching on</i>	26
6.2.4 <i>Shutdown</i>	26
6.3 Operating mode parameters	27
7 Rules for care, cleaning and disinfection	29
8 Marking	thirty
9 Packaging	34
10 Safety requirements and environmental protection	35
1 1 Transport and storage	36
1 2 Maintenance and repair	37
1 3 Disposal information	38
1 4 Manufacturer's warranties	39
1 5 Claim information	41
1 6 Electromagnetic compatibility	42
List of regulatory documents used by the manufacturer to confirm compliance	43

INTRODUCTION

The main purpose of the manual and this appendix to the instruction manual is to provide the user with instructions for operating the ESWT apparatus, as well as information regarding its repair and maintenance. The manual provides performance data for the ESWT device, its installation, operation and maintenance, as well as safety information. Before use, the user must carefully read the manual and this appendix to the instruction manual in order to operate the ESWT device correctly and without error, in order to achieve the specified safety and performance standards.

1 PURPOSE AND CONTRAINDICATIONS

Extracorporeal shock wave therapy (ESWT) machine model SALUS - Optimus - Pro **designed** for **high** intensity pulsed electromagnetic exposure and extracorporeal shock wave therapy. For more details, see "indications for use"

The ESWT device is suitable for professional use in hospitals and clinics by physiotherapists.

Scope : in hospitals and clinics .

Indications for use:

- stimulation of muscle fibers;
- improvement of local blood supply;

Neurology:

- Pain syndromes: dorsalgia , cervicalgia , trigemenalgia , tension headaches, compression syndrome of the greater occipital nerve, piriformis muscle syndrome, intercostal neuralgia, tunnel syndromes.

Polyneuropathy

- Peripheral paresis: paresis of the facial nerve, paresis of the peroneal nerve, paresis of the radial nerve, paresis of the median nerve, paresis of the ulnar nerve, paresis of the recurrent nerve.

Plexopathies

Intervertebral hernia

Paresis central

- spinal injury, consequences of a stroke of the brain and spinal cord, consequences of operations on the brain and spinal cord.

Dystonia

hypothesia

Rehabilitation after neurosurgical operations

Traumatology:

Joint diseases: arthrosis, arthritis

Rehabilitation: after arthroscopic operations, after injuries

Pain syndromes: impingement syndrome, myofascial syndrome, epicondylitis , plantar fasciitis , tibial tract syndrome , sacroiliac joint syndrome, trochanteritis , temporomandibular joint pathology, Achilles tendon pathology, ligamentitis

Movement disorders: flat feet, scoliosis, dystonia, inactivity atrophy syndrome

Lymphedema

Urology and gynecology:

- Overactive bladder
- Erectile dysfunction
- premature ejaculation
- Syndrome of acute and chronic pelvic pain
- Rehabilitation after prostate surgery
- Hypotension of the bladder
- Detrusor sphincter dyssynergy
- Hyposthesia of the glans penis
- Pudendal nerve pathology
- Stress urinary incontinence

Gastroenterology:

- Irritable Bowel Syndrome
- Dysphagia

Cosmetology:

Increased muscle tone: abdomen, gluteal region, shoulder, face

Decongestive therapy

Maxillofacial Surgery:

- Rehabilitation after CLS operations;
- shoulder syndrome;
- humeroscapular periarthritis ;
- physiotherapy for fractures ;
- chronic proximal plantar fasciitis in patients aged 18 years or older with symptoms for 6 months or more;
- plantar fasciitis ;
- traction degeneration of the plantar fascial band at the beginning of the medial tubercle of the calcaneus;
- repeated muscle contraction by passing electric currents through the electrodes without contact with the affected area of the body;
- removal of pain sensations (myofascial pain, thick fiber pain , fine fiber pain , neuropathic pain , phantom pain, fibromyogia);
- treatment of osteoporosis;
- treatment of fibromyalgia , arthritis, chronic fatigue.
- providing a completely non-invasive high-intensity pulsed electromagnetic effect on the musculature of the pelvic floor in order to rehabilitate weak pelvic muscles

and restore neuromuscular control in the treatment of urinary incontinence in humans.

Contraindications for use:

Contraindications for use in extracorporeal shock wave therapy are as follows:

- should not be used in conjunction with other electronic medical devices;
- **do not perform procedures in the area of the endoprosthesis;**
- The following patients must be approved by a doctor for the procedure:
 - patients with high body temperature (more than 37 °C), the elderly (over 85 years old) and children (under 3 years old);
 - person with implants (eg, artificial pacemaker, implantable infusion port, hearing aid).
- The following patients should not be treated with the product if:
 - blood disorders, clotting problems, or use of anticoagulants;
 - blood-thinning drugs (warfarin);
 - thrombosis;
 - oncological diseases;
 - with diabetes mellitus, diabetic neuropathy ;
 - pregnant women;
 - acute inflammation in the target area;
 - cartilage growth in children;
 - contact infectious patients (to avoid transmission of infection);
 - corticoid therapy for 6 weeks.


Attention! Shock waves should not be applied over air-filled tissues (lungs), nor should they be used near major nerves, vessels, vertebrae, or on the head.

Contraindications for use in high-intensity pulsed electromagnetic exposure are as follows:

- do not use this equipment in parallel with other electronic medical equipment;
- care should be taken that the magnetic effect does not penetrate into the region of the heart;
- Generally, the following categories of patients cannot be treated with this equipment (approval from the attending physician is required prior to any treatment with this equipment):
 - patients with high body temperature (more than 37 °C), pregnant women, the elderly (over 85 years old) and children (under 3 years old);
 - a patient with a history or diagnosis of epilepsy or convulsive disability;
 - a patient with suspected epilepsy based on electroencephalography;
 - patients with traumatic brain injury and spinal cord injury within a month from the moment of injury;
 - **patients with implantable infusion ports or hearing aids;**
 - patient with cranial implants.
- use during magnetic resonance imaging.

Possible complications:

If any abnormalities are found, such as skin redness, pain, hematoma with petechiae, hematoma, bleeding, swelling, you should immediately stop using the device and consult a doctor.

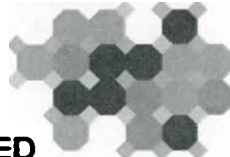
**REMED**

Possible side effects:

Side effects of extracorporeal shock wave therapy:

- edema;
- redness;
- hematoma;
- pain.

With pulsed electromagnetic exposure of high intensity, subject to the conditions and rules of transportation, storage and use, side effects have not been identified.

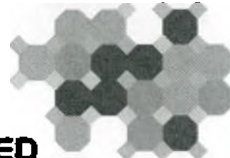REMED

2 GENERAL DESCRIPTION OF THE ESWT APPARATUS

The ESWT device consists of a body, a power cord, a shock wave therapy gun, an electromagnetic exposure transducer, and applicators for shock wave therapy.

The body of the ESWT device is equipped with a switch with a key, an emergency stop button, a touch LCD display (LCD -8 ") with a joystick , a compartment for storing applicators, rollers for movement with and without a brake, an electromagnetic impact converter, a holder lever assembly, a converter holder, connectors for connecting the shock wave gun and the power cord, as well as pre-installed software (SW).

The display of the ESWT device is made built into the case and is located on the top panel.

The shock wave gun consists of a body, a striker, a guide and a rubber insert, and is also equipped with a cable for connection to the body of the ESWT machine.

Shock wave applicators are available in the following sizes: 6mm, 10mm, 15mm, 20mm and 35mm.

3 DESIGN FEATURES AND APPEARANCE ESWT DEVICE

DESIGN FEATURES AND APPEARANCE OF THE ESWT DEVICE

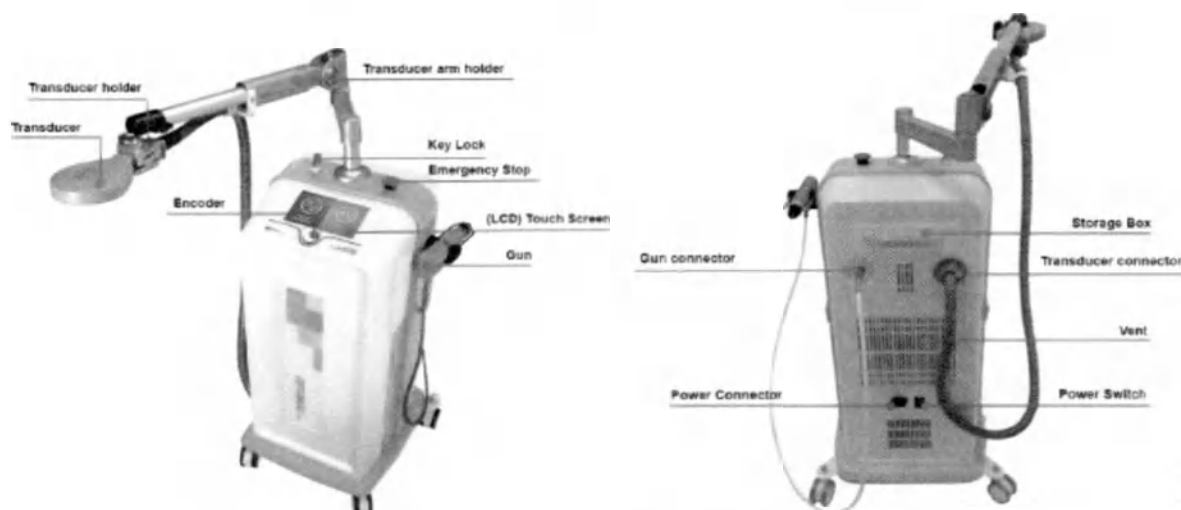
Figure 1. Appearance of the device in front (left) and rear (right)

- Transducer – electromagnetic impact transducer;
- transducer holder - the holder of the electromagnetic impact transducer;
- transducer arm holder - the lever of the holder of the electromagnetic impact transducer assembly;
- key lock - key switch;
- Emergency Stop – emergency stop button;
- Touch Screen LCD display (LCD-8");
- Encoder - joystick;
- Gun - gun shock wave therapy;
- Storage box - storage compartment;
- Gun connector - shock wave therapy gun connector;
- transducer connector – connector for connecting the electromagnetic impact transducer;
- Vent - ventilation holes;
- power connector - power cord connector with a holder bracket;

Figure 2. Power connector - power cord connector with holder bracket
power Switch - power switch.

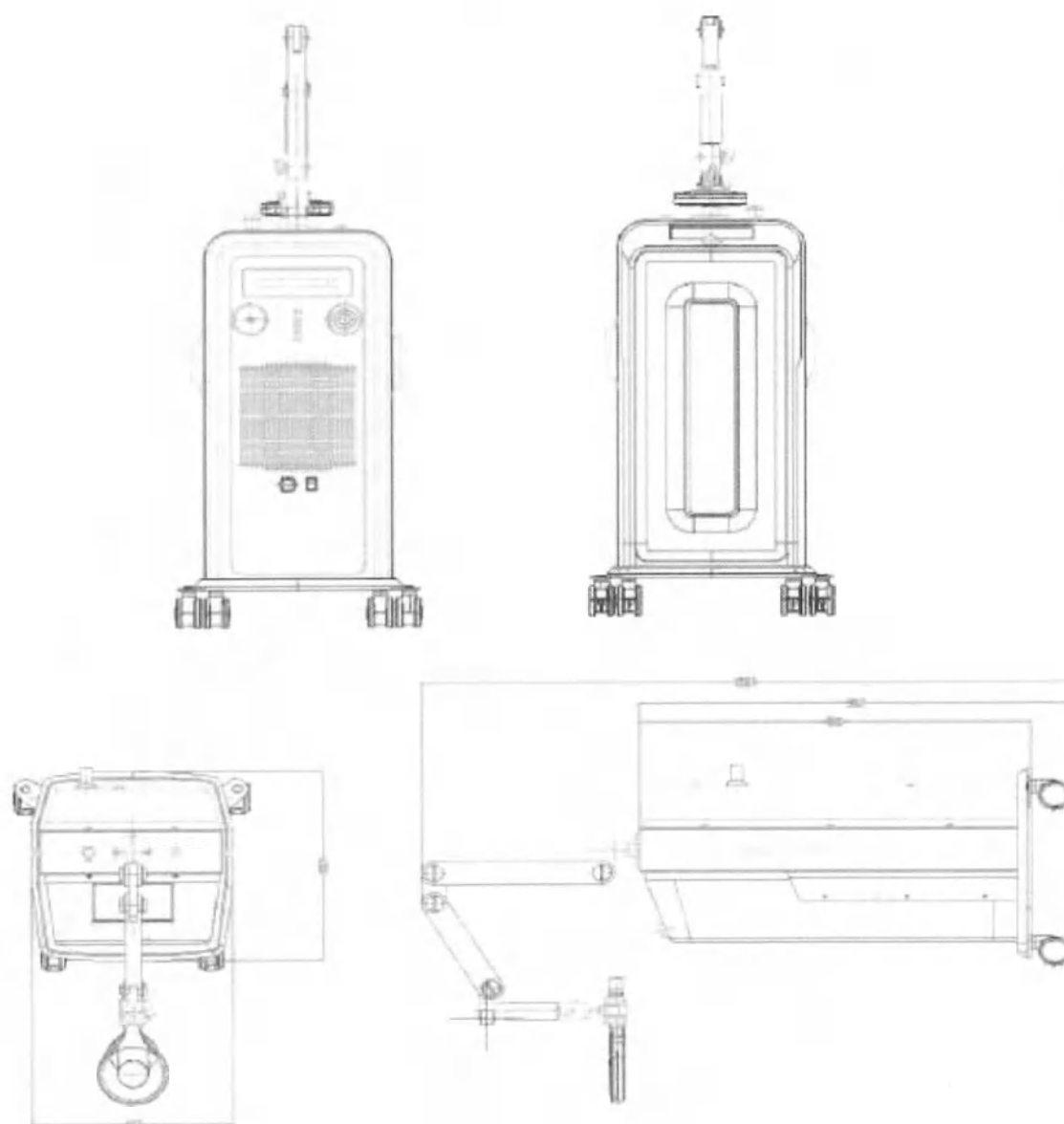
Figure 3. Dimensions of the ESWT apparatus

The body of the device is equipped with rollers with and without brakes for movement:

Figure 4. Roller with brake open (left) and closed (right)

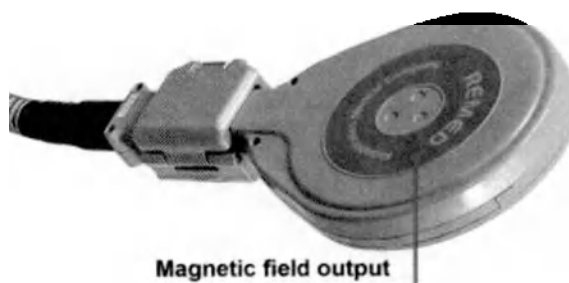
Figure 5. External view of the electromagnetic influence transducer
Magnetic field output is the surface of electromagnetic influence.

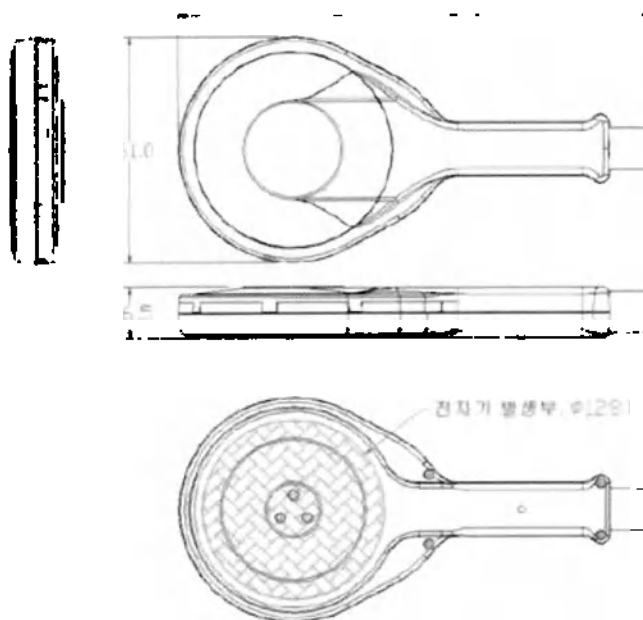
Figure 6. EMI transducer dimensions

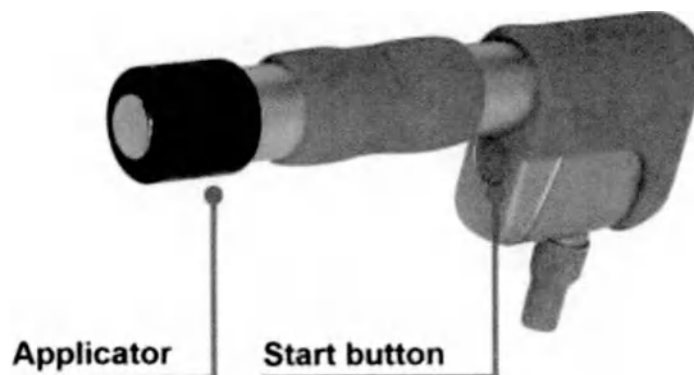
Figure 7. External view of the shock wave therapy gun
Applicator - applicator;
start button - start button.

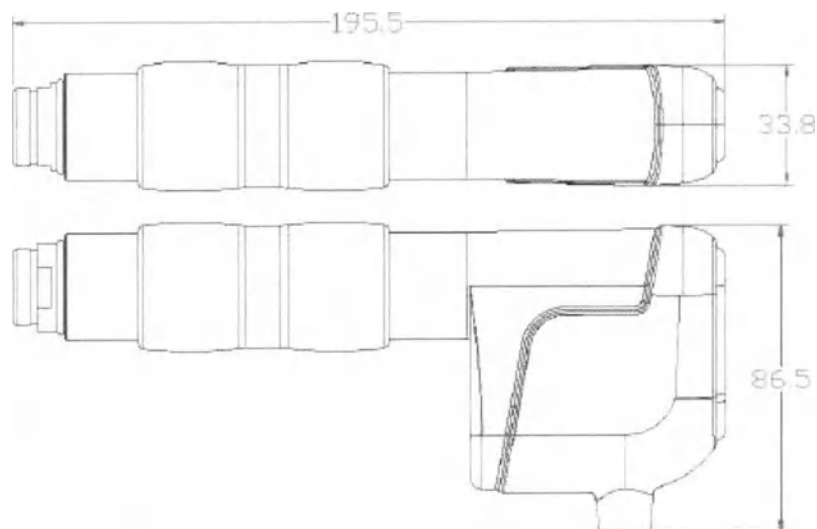
Figure 8. Shock Wave Gun Dimensions

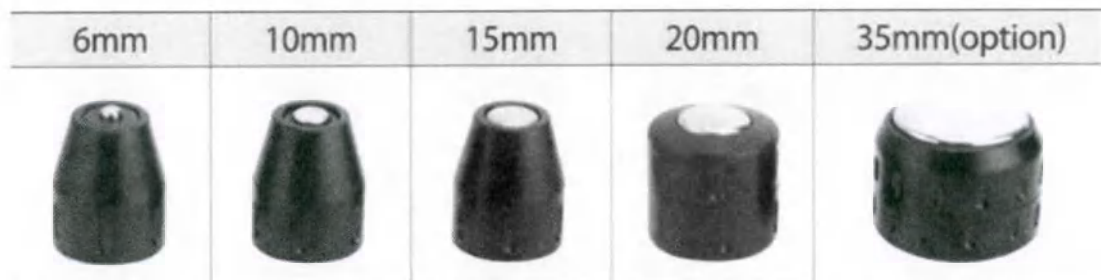
Figure 9. The appearance of the applicators and their size

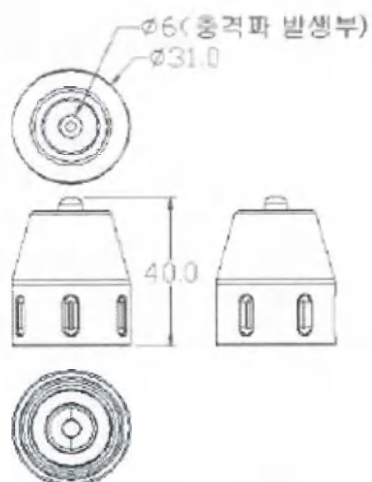
Figure 10. Applicator 6 mm

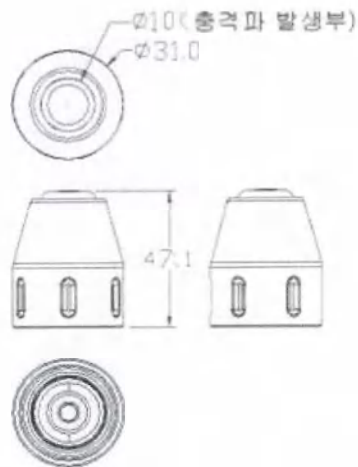
Figure 11 Applicator 10 mm

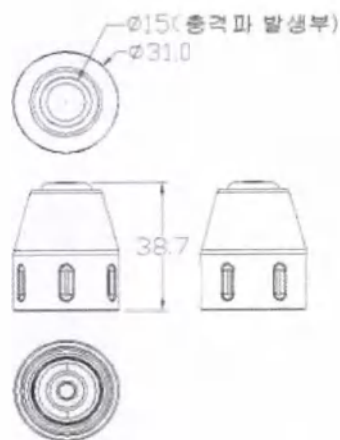
Figure 12. Applicator 15 mm

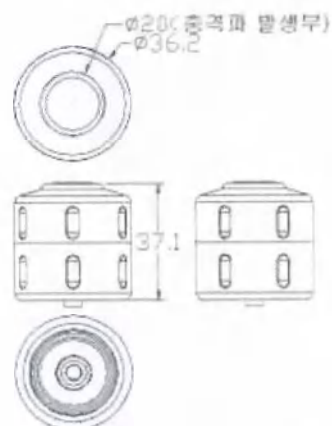
Figure 13 Applicator 20 mm

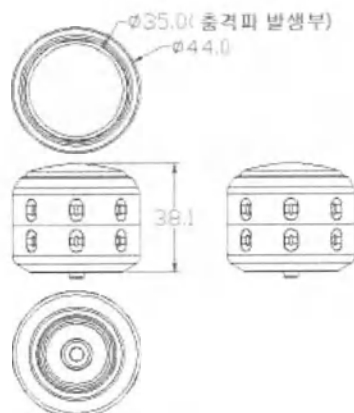
Figure 14 Applicator 35 mm

The ESWT device has a pre-installed proprietary software (hereinafter referred to as software). After turning on the product, the software starts automatically.

Table 1: Software details

Name of the software	Software identification name	Version number (identification number) of the software	Digital software identifier (checksum of the executable code)	Software Digital ID Calculation Algorithm	Date of first software release
Optimus-Pro F/W (firmware software)	-	1.00.39 _	-	-	09/19/2019
Optimus-Pro S/W (software)	-	2.01.03	-	-	09/19/2019

Note: Firmware and software perform different functions together in one product. Firmware is designed to control the equipment - controls the output according to the input from the software), the software is designed to control only what is shown on the display: therapy time, buttons.

Warning

- Software updates should only be carried out by an authorized representative. The software is encrypted and cannot be modified by yourself.

REMED

CONTENTS OF DELIVERY:

- 1) Apparatus ESWT , model SALUS-Optimus-Pro - 1 pc ;
- 2) Shock wave therapy gun with cable - 1 pc;
- 3) Power cord - 1 pc;
- 4) Applicator 6 mm - 1 pc;
- 5) Applicator 10 mm - 1 pc;
- 6) Applicator 15 mm - 1 pc;
- 7) Applicator 20 mm - 1 pc;
- 8) Applicator 35 mm - 1 pc;
- 9) Operation manual - 1 pc.

4 TECHNICAL DATA

The operational and technical characteristics of the ESWT device correspond to tables 2-3.

Table 2

Parameter name	Value, norm
1	2
Overall dimensions of the ESWT apparatus (H×L×W), mm	1960×515×580
ESWT body dimensions (H×L×W), mm	993×467×501
Dimensions of the compartment for storing applicators (H×L×W), mm	70×242×72
Roller diameter with brake, mm	75
Number of rollers with brake, pcs	2
Roller diameter without brake, mm	75
Number of rollers without brake, pcs	2
Mass of the body of the ESWT apparatus, kg	76
Shock wave therapy gun dimensions (H×L×W), mm	33.8×195.5×86.5
Shock wave therapy gun cable length, mm	1935
Shock wave therapy gun holder dimensions (H×L×W), mm	46×77×124
Number of shock wave therapy guns, pcs	1
Weight of shock wave therapy gun, kg	0.6
Dimensions of electromagnetic impact transducer (H×L×W), mm	36.5×330×161
Length of a hose of the electromagnetic influence converter, mm	2000
Number of electromagnetic impact transducers, pcs	1
Mass of electromagnetic impact converter, kg	3
Number of operating modes, pcs	2 (shock wave therapy and electromagnetic influence)
Time of readiness for work after inclusion, min., no more	1
Maximum time of continuous operation, min., not less than	Long-term with SWT Prolonged with EMW
Power cord length, mm	1500
Power cord weight, g	130
Rated supply voltage, V	220
Power consumption, kVA, no more	1.5
Nominal value of the mains frequency, Hz	50
Display type	built-in
Product weight without packaging / with packaging, kg	76/91
Joystick dimensions (Diameter×Height), mm	34×14
Emergency stop button dimensions (Diameter×Height), mm	40×28
Key switch base dimensions (Diameter×Height), mm	30×11

Key switch dimensions (Diameter×Height), mm	21×14
<i>Display specifications:</i>	
Display type	color TFT LCD (Touch Screen)
Dimensions (L×W), mm	180×103
<i>Characteristics of the holder lever with the holder of the electromagnetic impact transducer assembly:</i>	
Overall dimensions (H×L×W), mm	142×980×97
Holder lever base length, mm	295
Angle of rotation of the base of the holder lever around the vertical axis of the apparatus body, deg	360
Bracket dimensions for attaching the lever to the base of the lever (H×L×W), mm	75×105×71
Angle of rotation of the bracket with the lever around the vertical axis of the base of the holder lever, deg	360
Lever length, mm	553
Angle of rotation of the lever in the bracket, degrees	135
Bracket dimensions for attaching the transducer hose to the arm (H×L×W), mm	20×78×97
Bracket length for fixing the transducer holder to the arm, mm	78
Bracket diameter for fixing the transducer holder to the lever, mm	44
Dimensions of the screw handle at the bracket for attaching the converter holder to the bracket arm (H×L×W), mm	20×50×20
Transducer holder overall dimensions (H×L×W), mm	97×76×72
Internal dimensions of the transducer holder (H×L×W), mm	36×40.5×72
Dimensions of transducer holder axis (Diameter×Length), mm	12×12
Transducer holder base diameter, mm	31
Angle of rotation of the transducer holder around the axis of the lever, deg	360
Angle of rotation of the transducer holder around the axis perpendicular to the lever, deg	180

Note - Tolerances for all characteristics - $\pm 2.5\%$, unless otherwise indicated.

T a b l e 3 Characteristics of applicators

Applicator size	Dimensions (heightxdiameter), mm	Impact surface diameter, mm	Weight, g
1	2	3	4
Applicator 6 mm	40.0x31.0	6	31
Applicator 10 mm	47.1x31.0	10	42.5
Applicator 15 mm	38.7x31.0	15	57
Applicator 20 mm	37.1x36.2	20	68
Applicator 35 mm	38.1x44.0	35	100.5

N o t e - Tolerances for all characteristics - $\pm 2.5\%$, unless otherwise indicated.

The ESWT apparatus remains operational at temperatures from $+10^{\circ}\text{C}$ to $+30^{\circ}\text{C}$, relative air humidity from 30 to 85% and atmospheric pressure from 70 to 106 kPa.

The operating mode is continuous and depends on the selected parameters.

Protection class against electric shock: I , applied part Type BF .

Protection against ingress of water and dust IP X0.

Service life - 3 years.

5 MATERIALS USED IN MANUFACTURING

ESWT devices are made of the following materials:


- ABS plastic PA 200 – for body and body elements, shock wave therapy gun, applicators;
- Stainless steel 420 SS - for the body and body elements;
- XINJIANG grade SG-3 polyvinyl chloride for power cord, gun cable and transducer hose;
- Rubber grade IR -307 - for inserting the shock wave therapy gun and applicators;
- Stainless steel grade AMS 5355 by Wire Company , Inc. _ - for applicators;
- Polycarbonate brand Lupoy ® GN1004FA manufactured by LG Chem ltd . – for the converter of electromagnetic influence.

ESWT devices do not contain drugs, components of animal or human origin.

In the manufacture of ESWT devices, urethane dyes are used: WHITE, BLACK, LIGHT VIOLET (PANTONE 270C), COOL GRAY (PANTONE 10C).

All color options are provided by mixing basic dyes.

Type and duration of contact with the human body: Short-term contact (no more than 24 hours) with the surface of intact skin.

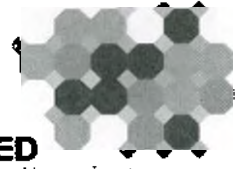REMED

6 USING THE ESWT MACHINE

6.1 Unpacking and installation

Open the package in accordance with the symbols indicated on it and carefully remove the ESWT device and its components:

- Check the ESWT apparatus and its components for any mechanical damage;
- In case of temperature difference, it is necessary that the ESWT apparatus be at room temperature for at least 8 hours;
- If any problems are found, contact the manufacturer's authorized representative immediately.

⚠ Warning

- **Save packing materials for future transport and storage;**
 - **Disposal of packaging materials must be carried out in accordance with your local regulations.**
-

6.2 Preparing the ESWT apparatus for operation

6.2.1 Installation method

When moving the ESWT apparatus, rollers on the body are used;

To move, set the braked castors to the open position (see Figure 4);

After moving, set the brake rollers to the closed position (see Figure 4);

Place the ESWT machine on a flat and level surface;

Install the ESWT in a place with suitable ambient temperature and humidity, and do not install in a place where it will be exposed to dust or flammable materials. Operate at temperatures from +10°C to +30°C and relative air humidity from 30 to 85%;

Connect the power cord to the body of the ESWT (see Figure 2);

Connect the power cord to a grounded socket;

Make sure the power cord is connected;

Before turning on the equipment, make sure that the socket and power cable are grounded;

Be careful not to damage the equipment with strong impacts;

Plug in the shock wave therapy gun;

Fix the electromagnetic influence transducer in the holder. When disconnecting the transducer from the holder, hold it with your hand so that when the holder lever is lifted with a gas shock absorber, it does not fall;

Turn on the power switch located on the back of the ESWT (see Figure 1);

When the power is turned on, the manufacturer's logo will appear on the display of the ESWT device for a moment and then the "Main Screen" appears, on which you can select the operating mode;

The patient is seated or placed in front of the apparatus;

Connect the appropriate applicator to the shock wave therapy gun;

The transducer or applicator is placed over the desired area of the body depending on the pain;

Choose the necessary program or set the necessary parameters for the procedure;

During the procedure, the display shows data on the parameters of shock wave therapy or electromagnetic exposure, as well as warnings;

The duration of the procedure depends on the selected parameters.

Power connection

Please use the original power cord.

Figure 16. Power cord

Please make sure the AC power supply meets the following specification: AC 230(V), 50(Hz).

Connect the power cord to the power cord connector on the back of the ESWT. Connect the opposite end of the power cord to an AC outlet suitable for use with medical equipment, the outlet must be grounded. If you are in doubt about the availability of grounding, contact the appropriate personnel in your healthcare facility.

Note

- **Lack of grounding leads to a breakdown of static electricity and the failure of the ESWT apparatus.**

6.2.3 Switching on

After installing the ESWT device, please follow the procedures below to turn it on:

- Before using the ESWT apparatus, carry out an appropriate safety check;
- Press the power switch on the back of the ESWT machine;
- Turn the key switch to the "ON" position on the top panel of the ESWT machine;
- The boot image will appear on the screen;
- After 12-15 seconds, the system will display "Home Screen";
- Now you can operate the ESWT device using the touch screen and control joystick.

6.2.4 Shutdown

To turn off the ESWT device, please follow the procedures below:

- Turn the key switch to the "OFF" position on the top panel of the ESWT;
- Press the power switch on the back of the ESWT.

6.3 Operating mode parameters

The parameters of the operating modes correspond to tables 4-5.

Table 4

Parameter name	Value, norm
1	2
Shock Wave Therapy Options:	
Operating modes	Manual with the possibility of entering into the memory of the device 20 programs with adjusted parameters of shock wave therapy (from M 1 to M 20) Adjustable parameters: number of strokes, frequency, pressure
Frequency range, Hz with resolution, Hz and tolerance, %	1-21 1 ± 20
Pressure range, bar with resolution, bar and tolerance, %	1-6 0 ± 20
Range of number of beats, beats with discreteness, strokes	100-5000 100
Standard sizes of used applicators	6, 10, 15, 20, 35
Parameters of electromagnetic influence:	
Operating modes	Automatic with the ability to select one of 20 programs with already set electromagnetic exposure parameters, and with the ability to adjust the time and magnetic field strength (from A1 to A20) (see Table 5) Manual with the possibility of entering into the memory of the device 20 programs with adjusted parameters of electromagnetic influence (from M 1 to M 20). Adjustable parameters: frequency 1, frequency 2, exposure time and break time. User with the possibility of entering into the memory of the device 20 programs with adjusted parameters of electromagnetic influence (from

	U 1 to U 20). Adjustable parameters: frequency, exposure time, break time, magnetic field strength and repetition cycles.
Frequency range, Hz	1-100
with resolution, Hz	1
and tolerance, %	±20
Exposure time range, s	1-4
Time range of breaks between impacts, s	0-10
Procedure time range, min	1-60
with resolution, min	1
Range of repetition cycles in the procedure, cycles	1-100

Table 5

Purpose		Program	Time (min)	frequency Hz
acute pain	Stretching , irritation	A1, A2	10 - 15	30-50
	Ankle pain	A3	10 - 15	5-7
	Shoulder pain	A4	10 - 15	3-5
	Pelvic/thigh pain	A5	10 - 15	30-50
	humeroscapular periarthrosis	A6	10 - 15	5-10
	Backache	A7, A8	10 - 15	30-50
	Arthritis, knee	A9	10 - 15	5-7
	Elimination of pain	A10	10 - 15	5-7
chronic pain	Stretching , irritation	A11, A12	10 - 15	10-20
	Ankle pain	A13	10 - 15	5-7
	Shoulder pain	A14	10 - 15	1-5
	Pelvic/thigh pain	A15	10 - 15	10-20
	humeroscapular periarthrosis	A16	10 - 15	5-7
	Backache	A17, A18	10 - 15	10-20
	Arthritis, knee	A19	10 - 15	3-5
	Elimination of pain	A20	10 - 15	3-5

7 RULES FOR CARE, CLEANING AND DISINFECTION

⚠ Warning

- Before cleaning the gun/transducer, turn off the ESWT and disconnect the gun cable from the connector;
- Preventive checks and cleaning of the applicators, gun, transducer and external surfaces of the ESWT apparatus body can be carried out by users;
- Cleaning of the inside of the ESWT body is only permitted by personnel authorized by the manufacturer.

Your ESWT device should be cleaned regularly. The outside of the instrument can be cleaned using a clean and soft cloth moistened with 70% isopropyl alcohol.

Disinfection can cause equipment damage, only do it when included in the maintenance plan and clean the ESWT before disinfecting.

Disinfectant - 70% isopropyl alcohol.

⚠ Warning

- NEVER immerse the ESWT device and its parts in water or any solution, pour water or any solution on the ESWT device and its parts;
- NEVER allow liquid to enter the housing, switches, connectors, or any vents on the ESWT;
- ALWAYS wipe off excess liquid remaining on the surface of the ESWT apparatus and its parts with a cloth;
- NEVER use abrasive, corrosive cleaners or cleaners containing acetone;
- NEVER disinfect the ESWT device and its components using high pressure and high temperature.

Periodically clean the exterior of the ESWT machine, gun, and touch screen with a soft cloth dampened with 70% isopropyl alcohol.

After each shock wave treatment, clean the applicator with 70% isopropyl alcohol and a smooth cloth. Regular cleaning of the applicators should be carried out once a week, as the gel used during the procedures may stick to the applicator after prolonged use and may affect the generation of the shock wave.

Use running water to clean the applicators. If the gel does not wash off well, it is effective to submerge the applicator in water for 2-3 hours and then rinse. After washing, remove water with a soft dry cloth and dry at room temperature for more than 1 hour.



To ensure safe use, the internal cleaning of the ESWT housing must be carried out once a year by personnel authorized by the manufacturer.

8 MARKING


8.1 Marking applied to the label affixed to the body of the ESWT apparatus:

- ESWT device model;
- manufacturer's logo;
- serial number;
- date of manufacture;
- Rated voltage;
- alternating current;
- rated frequency;
- maximum power consumption;
- authorized representative in the European Community;
- working part type BF
- general warning sign;
- refer to the instruction manual;
- do not dispose of with household waste;
- name and address of the manufacturer;
- do not push;
- do not sit;
- CE marking.

8.2 Safety labels affixed to the surface of the apparatus:

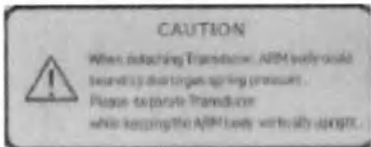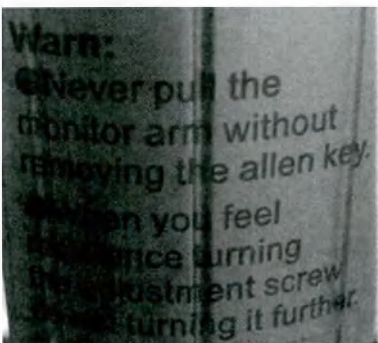Label	Place of gluing	Description
	Upper side surface of transducer holder arm	Warning sign related to the transducer holder lever: CAREFULLY When the transducer is disconnected, the holder arm may move up due to the pressure of the gas buffer. Please keep the holder lever vertical when disconnecting the transducer.
	The upper surface of the lever of the converter holder at the bracket for attaching the lever to the base of the lever	Warning sign related to the transducer holder lever: Warning: • Never pull the holder lever without removing the hex wrench. • When you feel resistance when turning the adjusting screw, do not turn it any further.


8.3 Marking applied to the packaging of the device:

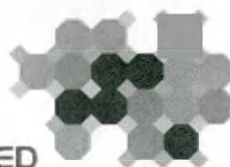Symbol	Meaning	Description
	Do not stack.	Indicates that items cannot be stacked vertically, either due to the nature of the shipping packaging or the nature of the items themselves.
	Keep out of the rain	Indicates that the transport packaging should be stored away from rain and in dry conditions.
	Fragile	To indicate that the contents of the shipping package are fragile and should be handled with care.
	This side up	To indicate the correct vertical position of the transport package.
	Recycling	To specify the location of a trash can or container

8.4 Decoding of all marking symbols:

Symbol	Meaning	Description
	Warning	Warning of danger to the life of the operator if the precautions are not observed
	Attention	Notice of precautions to be taken to avoid injury or damage
	Don't push	It is forbidden to push
	Don't sit	It is forbidden to sit on the surface
	Type BF	Working part type BF
	Work instructions	You must follow the operating instructions
	Operating instructions	You must follow the operating instructions
	Serial number	Product serial number
	Authorized Representative in the European Community	Authorized Representative in the European Community

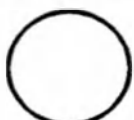Symbol	Meaning	Description
	vendor code	Product number
	Manufacturer	Manufacturer
	date of manufacture	date of manufacture
	Do not stack	Note on the prohibition of stacking products in packaging
	Keep away from moisture	Instruction on the need to store the transport packaging away from rain and in dry conditions
	fragile	Fragile. Carefully
	Top	Specifying the Correct Vertical Position of the Shipping Package
	Reuse	Indication of the need for disposal at the appropriate collection points for waste to be recovered or recycled
	Off (power)	Designation of the position of the mains switch corresponding to disconnection from the mains
	Included (food)	Designation of the position of the mains switch corresponding to the inclusion in the network
	Alternating current	Indication on the product label about the need to connect the product through the appropriate connector only to AC sources
	Do not dispose of with household waste	Indication of the need for disposal at separate collection points for recovery and recycling
	Temperature Range	Indication of the maximum and minimum temperatures for storage, transport and operation

Symbol	Meaning	Description
	Humidity range	Specify the permissible upper and lower limits of relative humidity during transportation and storage
	Atmospheric pressure range	Indication of permissible upper and lower limits of atmospheric pressure during transportation and storage
	dangerous voltage	Caution, dangerous voltage
	CE marking	CE marking

9 PACKAGING

The stimulator, packed in a plastic bag, is placed in shock-absorbing foam inserts. On top of the shock-absorbing insert, in a special recess, a cardboard box is placed in which there is a power cord, bracket, instruction manual in Russian, keys. From above it is closed by a package from polyethylene. Then the outer carton is put on.

10 SAFETY REQUIREMENTS AND ENVIRONMENTAL PROTECTION

10.1 Electrical safety precautions:

- The rated voltage of the power supply for the ESWT apparatus is 220 V, alternating current;
- Make sure that all connecting parts (power cord or component cables) are properly connected to the body of the ESWT;
- Make sure the ESWT apparatus is fully grounded;
- Repair, expansion and installation of the ESWT apparatus must not be carried out by anyone other than specialized personnel authorized by the manufacturer. Arbitrary disassembly / assembly of the ESWT device by the user is strictly prohibited;
- Before connecting other equipment not listed in the instruction manual, be sure to consult the manufacturer or an authorized representative of the manufacturer;
- To avoid electrical interference during use, the ESWT device should be installed at a considerable distance from the power source, any generator, x-ray equipment, transmitter, mobile electrical wire, etc.;
- Be careful of possible electromagnetic or other interference between the ESWT machine and other devices;
- As such, an independent power circuit is required and power sharing with other electronic devices is not recommended;
- Do not place the ESWT in a location where it is difficult to disconnect the power plug or other detachable plug.

For safety, the ESWT apparatus must comply with the requirements of IEC 60601-1.

The ESWT model SALUS - Optimus - Pro must additionally comply with the requirements of IEC 60601-2-10.

For electromagnetic compatibility, the ESWT apparatus must comply with the requirements of IEC 60601-1-2 (class A, group 1).

Serviceability in accordance with IEC 60601-1-6 and IEC 62366.

Software class according to IEC 62304 - B.

10.2 Environmental protection

The materials used in the manufacture of the ESWT device should not pose a danger to human life, health and the environment both during operation and after the end of its service life.

The main type of possible hazardous impact on the environment is pollution of the atmospheric air of populated areas, soils and waters as a result of:

- emergency leaks (placers) of production materials;
- unorganized waste disposal on the territory of the manufacturer or outside it;
- random dumping of waste in places not intended for this purpose.

11 TRANSPORT AND STORAGE

11.1 Transport conditions

Transportation of the ESWT device can be carried out by all types of covered vehicles and the rules for the carriage of goods in force for each type of transport.

Transportation of the ESWT device by sea should be carried out in accordance with the "Safety Rules for the Sea Transportation of General Cargoes". Type of shipment - containers and small shipment.

The ESWT device, during transportation, must be resistant to climatic influences at an ambient temperature of -10°C to $+60^{\circ}\text{C}$, relative humidity of air not more than 80% and atmospheric pressure from 70 to 106 kPa.

The ESWT device, during transportation, must be resistant to mechanical stress, with the following parameters:

- vibration loads - frequency range 10-55 Hz with a displacement amplitude of 0.15 mm;

During transportation it is prohibited:


- push the package with the ESWT apparatus;
- sit on the package with the ESWT apparatus;
- stack packages with ESWT devices.

11.2 Storage conditions

The ESWT device, when stored, must be resistant to climatic influences at an ambient temperature of -10°C to $+60^{\circ}\text{C}$, relative air humidity of not more than 80% and atmospheric pressure of 70 to 106 kPa.

Storage conditions:

- Keep away from water;
- Keep away from direct sunlight;
- Do not store near heaters;
- Do not store in places where dust accumulates;
- Avoid locations subject to strong shock or vibration, chemicals or explosive gases.

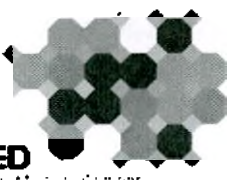REMED

12 MAINTENANCE AND REPAIR

⚠ Warning

- The user should routinely inspect the ESWT device, cables, gun, transducer, and applicators for damage that could affect patient safety or performance. The recommended inspection interval is one week or less. Repair damage before using the instrument;
- Periodic safety checks must be carried out by personnel authorized by the manufacturer. The recommended interval is one year, or you must follow the hospital's established intervals and procedures.

12.1 Routine inspection of the ESWT apparatus:

- The outer sheath of the ESWT power cord, gun cable/transducer hose, etc. must not peel off, the internal wires must not be exposed or damaged by external impacts;
- There should be no evidence of oil leakage from the converter;
- Clean the outside of the ESWT apparatus so that there are no foreign objects;
- The joystick for controlling the ESWT apparatus, etc., must not tremble;
- The various parts attached to the ESWT apparatus must not shake.

12.2 After a long period of storage, be sure to check the ESWT device before using it.

12.3 Removing condensate from the air compressor

Condensation must be removed after every 1,500,000 discharges of use. Otherwise, it may cause a malfunction or an unpleasant odor. More frequent condensate draining may be required depending on the ambient humidity.

12.4 Replacing gun consumables

The hammer, guide and rubber insert of the gun are recommended to be changed after every 1,500,000 to 3,000,000 (3 bar, 15 Hz) strokes of use, but the replacement period may vary depending on operating conditions.

12.5 To ensure and guarantee safe use, be sure to check the ESWT apparatus, including internal components and output voltage, by personnel authorized by the manufacturer once a year.

⚠ Warning

Safety checks or maintenance that requires opening the casing of the ESWT device should only be carried out by an authorized representative of the manufacturer. Failure to do so may result in equipment failure and possible health risks .

Maintenance and repair of the ESWT device by the user is not allowed.

Maintenance and repair is carried out only by an authorized representative of the manufacturer.

It is forbidden to carry out repair work with the power supply switched on.

Infomed-Neuro Limited Liability Company (Infomed-Neuro LLC) at the address: Russia, 125445, Moscow, Leningradskoe shosse, bld. 69, building 1, room I, fl . 1 room 50; phone +7 495 645 47 00; email kitayeva @ im - neuro . ru .

13 DISPOSAL INFORMATION

Disposal of the ESWT device must be carried out in accordance with SanPin 2.1.3684-21 for waste class A, including applicators and the transducer after preliminary disinfection, belongs to class A - epidemiological safe waste.

The ESWT device must be disposed of if:

- end of service life;
- confirmation of facts and circumstances that pose a threat to the life and health of health workers and indicate failure to fulfill the intended purpose.

Recycling is also subject to all packaging, including transport. Paper, wood, polyethylene and plastic must be disposed of separately in accordance with the rules for the collection, accounting and disposal of packaging.

Do not disassemble the ESWT device except for the purpose of disposal!

Premature disassembly of the ESWT apparatus will lead to its destruction!

The ESWT device must be disconnected from the power supply before disassembly!

14 MANUFACTURER WARRANTY

The ESWT apparatus is guaranteed for a period of 12 months from the date of sale from the manufacturer's warehouse.

It is guaranteed that the ESWT apparatus is free from defects related to the production and quality of production materials. Any defective part will be replaced if the ESWT device has not been tampered with and the ESWT device has been properly used during this period. Make sure that any malfunction is not due to inadequate maintenance of the ESWT apparatus or failure to follow the instruction manual. In the event of a malfunction, subject to the operating requirements in accordance with the operating manual, the manufacturer performs maintenance or replacement of the ESWT apparatus.

The manufacturer or an authorized representative of the manufacturer in the territory of the Russian Federation shall not be liable for any personal / material damage caused by erroneous use, operation for purposes other than the intended purpose, and negligent maintenance of the product.

Please note that if the ESWT device is damaged using unauthorized cleaning and disinfection materials, no warranty is provided.

Authorized representative of the manufacturer in the territory of the Russian Federation: Limited Liability Company " Infomed-Neuro " (LLC " Infomed-Neiro ") at the address: Russia, 125445, Moscow, Leningradskoe shosse, 69, building 1, room. I, fl . 1 room 50; phone +7 495 645 47 00; email kitayeva @ im - neuro . ru .

WARRANTY CARD

Product category _____

Manufacturer _____

Version _____

Serial number/art. _____

Date of sale _____

Organization of the seller _____

Seller's signature _____

Stamp
about sale

The product has been tested and has no damage.

I have read and agree to the terms of use.

Buyer's signature _____

The warranty period is 12 months from the date of sale from the manufacturer's warehouse.

For warranty issues, please contact: Limited Liability Company " Infomed-Neuro " (LLC " Infomed-Neuro ") at the address: Russia, 125445, Moscow, Leningradskoe shosse, 69, building 1, pom. I, fl . 1 room 50; phone +7 495 645 47 00; email kitayeva @ im - neuro . ru .

Infomed-Neiro Limited Liability Company (Infomed-Neiro LLC) at the address: Russia, 125445, Moscow, Leningradskoe shosse, 69, building 1, pom. I, fl . 1 room 50; phone +7 495 645 47 00; email kitayeva @ im - neuro . ru .

ATTENTION! Make sure that the warranty card is filled out correctly. The serial number of the product / article must exactly match the number in the coupon. All columns of the coupon (with the exception of the columns for the guarantee center) must be filled in.

In case of inaccurate or incomplete filling of the coupon, the warranty for the product is considered lost.

The manufacturer's warranty is carried out in accordance with the Civil Code of the Russian Federation on a general basis.

The product must not bear traces of mechanical or thermal effects, as well as traces of exposure to various chemical agents. The product may only be operated in accordance with the instruction manual for the product.

WARRANTY STATEMENT

DATE	DESCRIPTION	WARRANTY CENTER

15 COMPLAINT DETAILS

In case of warranty and out-of-warranty repair of the product, contact: Limited Liability Company " Infomed-Neuro " (LLC " Infomed-Neuro ") at the address: Russia, 125445, Moscow, Leningradskoe shosse, 69, building 1, pom. 1, fl . 1 room 50; phone +7 495 645 47 00; email kitayeva @ im - neuro . ru .

16 ELECTROMAGNETIC COMPATIBILITY

1. EMC specifications

Definition: EMC (Electromagnetic Compatibility) is the ability of equipment to avoid electromagnetic interference from other devices without causing similar electromagnetic interference to other equipment.

Warning

Because this ESWT device generates, uses, and radiates radio frequency energy, other medical and non-medical devices and wireless communications may be affected. It is characterized by the creation of electromagnetic interference in the air or cable environment. This product is designed in full compliance with EMC standards and provides reasonable protection against such interference.

Warning

ESWT devices require special measures to ensure ELECTROMAGNETIC COMPATIBILITY and must be installed and put into service in accordance with the EMC-related information given in the instruction manual.

Note:

- Using devices such as cellular phones, transceivers, or remote control devices that emit radio waves near this product may interfere with the operation of this product;
- When installing, make sure that this equipment is located as far as possible from other electrical equipment;
- Use only cables provided or designed by our company;
- Do not attempt to modify this equipment. Modifications to the product may result in reduced EMC performance. Modifications include cable modifications, improper installation or location of the product, changes to its configuration or components, changes to established procedures for handling the device or accessories, etc.

2. EMC Corrective Measures

If this ESWT device is found to cause interference, the user is encouraged to correct the problem by one or more of the following methods:

- Keep other devices as far away from this ESWT as possible to reduce electromagnetic interference;
- To reduce electromagnetic interference, you can adjust the position of equipment and other devices or the angle between them;
- Electromagnetic interference can be reduced by changing the connection point of the ESWT power cord;
- You can also change the power supply path to other devices to reduce electromagnetic interference.

The manufacturer assumes no responsibility for interference caused by using a non-recommended connecting cable, or by modifying or altering the ESWT without authorization. Unauthorized modification or alteration of the ESWT device may void the user's authority to operate this ESWT device.

The list of regulatory documentation used by the manufacturer to confirm compliance:

EN 1041 , Information provided by the manufacturer of medical devices
ISO 13485 , Medical devices - Quality management systems - Requirements for regulatory purposes
ISO 14971 , Medical devices - Application of risk management to medical devices
ISO 15223-1 , Medical devices - Symbols to be used on medical device labels, markings and information to be provided - Part 1: General requirements
IEC 60601-1 , Medical electrical equipment - Part 1: General requirements for basic safety and essential performance
IEC 60601-1-2 , Medical electrical equipment - Part 1-2: General requirements for basic safety and essential performance - Related standards: Electromagnetic compatibility: Requirements and tests
IEC 60601-1-6 , Medical electrical equipment - Part 1-6: General requirements for basic safety and essential performance - Collateral standard: Ease of use
IEC 60601-2-10 , Medical electrical equipment - Part 2-10: Particular requirements for the safety of nerve and muscle stimulators
IEC 62304 , Medical device software - Software life cycle processes
IEC 62366 , Medical devices - Serviceability applications
MEDDEV 2.7/1 Revision 4 , CLINICAL EVALUATION: GUIDANCE FOR MANUFACTURERS AND NOTIFIED AUTHORITIES ACCORDING TO DIRECTIVES 93/42/ EEC and 90/385/ EEC
ISO 10993-1 , Biological evaluation of medical devices - Part 1: Evaluation and testing as part of the risk management process
ISO 10993-5 , Biological evaluation of medical devices - Part 5: Tests for cytotoxicity by the in method vitro
ISO 10993-10 , Biological evaluation of medical devices - Part 10: Skin irritation and sensitization tests
EN 980 - Symbols used on the labeling of medical devices

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