

Фотографии общего вида
«Пескоструйный аппарат DENTO-PREP Dental Microblaster с принадлежностями»

1. Пескоструйный аппарат DENTO-PREP Dental Microblaster



2. Распыляющая насадка/ Spray Tip



3. Разъем и шланг длиной 3 м / Quick Disconnect Coupling + 3 m hose



4. Разъем / Quick Disconnect Coupling




5. Шланг / Hose



6. Зарядная емкость для песка / Powder Jar




7. Крышка емкости для песка / Lid for Powder Jar

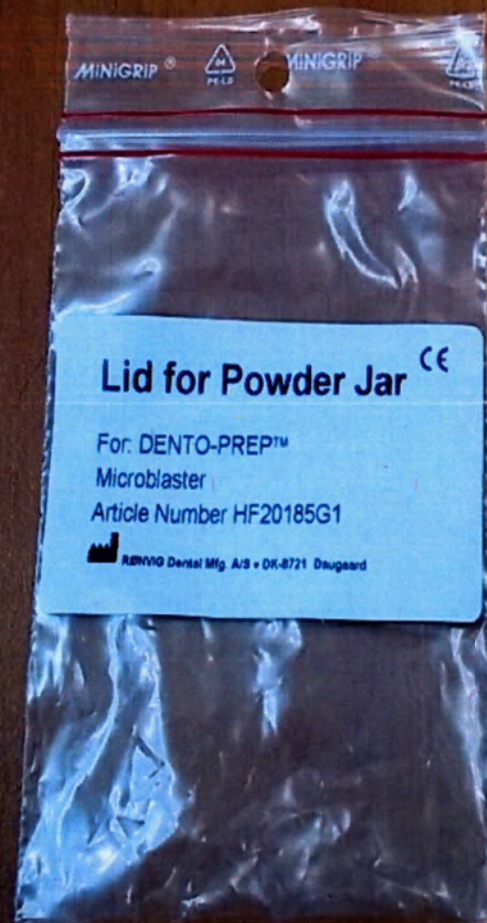

8. Разъем и фитинг для шланга диаметра 5 мм / Quick Disconnect Coupling and fitting for 5 mm hose




9. Разъем и фитинг для шланга диаметра 6 мм / Quick Disconnect Coupling and fitting for 6 mm hose

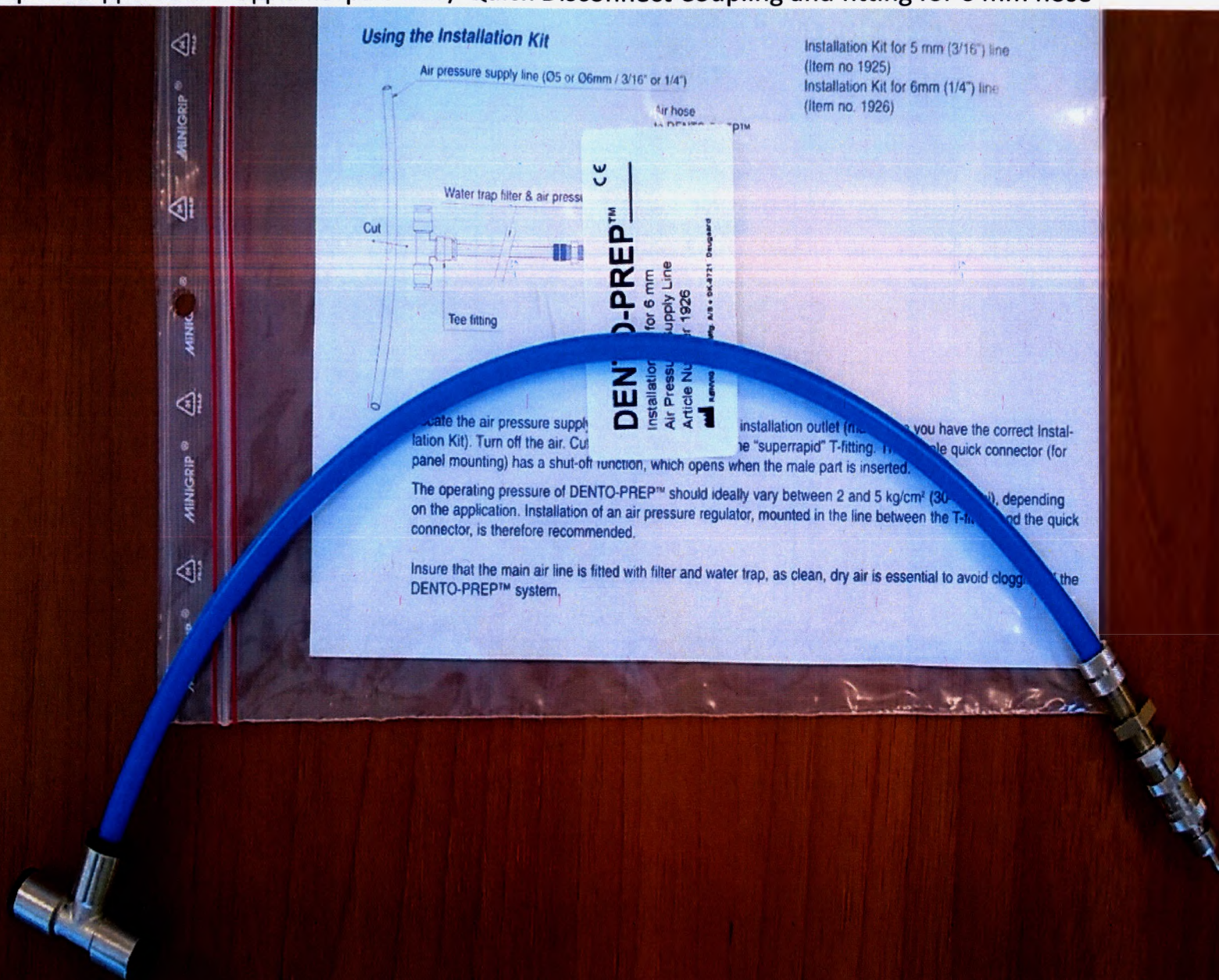

10. Разъем и фитинг для шланга диаметра 9 мм / Quick Disconnect Coupling and fitting for 9 mm hose

Using the Installation Kit

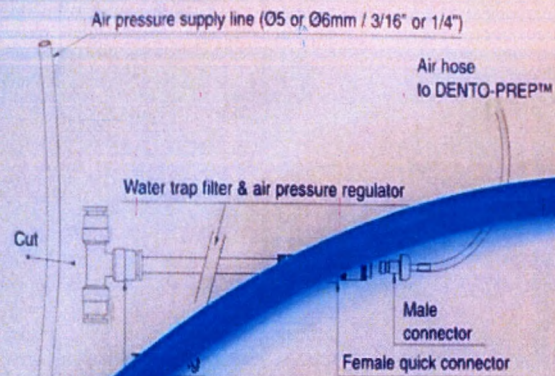

Installation Kit for 5 mm (3/16") line
(Item no 1925)

Installation Kit for 6mm (1/4") line
(Item no 1926)

DENTO-PREP™ CE
Installation Kit for 5 mm or 6 mm Air Pressure Supply Line
Article Number 1925 or 1926

Locate the air pressure supply line close to the desired installation outlet. (See the correct Installation Kit). Turn off the air. Cut the air hose and install the "superrapid" T-fitting. The female quick connector (for panel mounting) has a shut-off function, which opens when the male part is inserted.

The operating pressure of DENTO-PREP™ should ideally vary between 2 and 5 kg/cm² (30-120psi), depending on the application. Installation of an air pressure regulator, mounted in the line between the T-fitting and the quick connector, is therefore recommended.

Insure that the main air line is fitted with filter and water trap, as clean, dry air is essential to avoid clogging of the DENTO-PREP™ system.

11. Коннектор Midwest / Midwest Connector



12. Коннектор «S» для Sirona / Connector «S» for Sirona

Connector/S

For: DENTO-PREP™

Microblaster

Article Number 1935

 **Dewag** Dental Mfg. A/S • DK-8721 Degegaard

CE



13. Коннектор «W» для W&H / Connector «W» for W&H

DK-8721 Dagaard

Connector/W

CE

For DENTO-PREP™

Microblaster

Article Number 1936



DENTO-Prep Dental Mfg. A/S • DK-8721 Dagaard

CE



14. Коннектор «К» для Kavo / Connector «K» for Kavo

DK-8721 Daugaard

Connector/K

CE

For: DENTO-PREP™

Microblaster

Article Number 1937



DENTO Denta Mig A/S - DK-8721 Daugaard

CE



16. Коннектор «N» для NSK / Connector «N» for NSK





В настоящем документе
прошито и пронумеровано
16 (Шестнадцать) листов
Ген. директор ООО «АРКОМ»

Селезнев Р.Ю.