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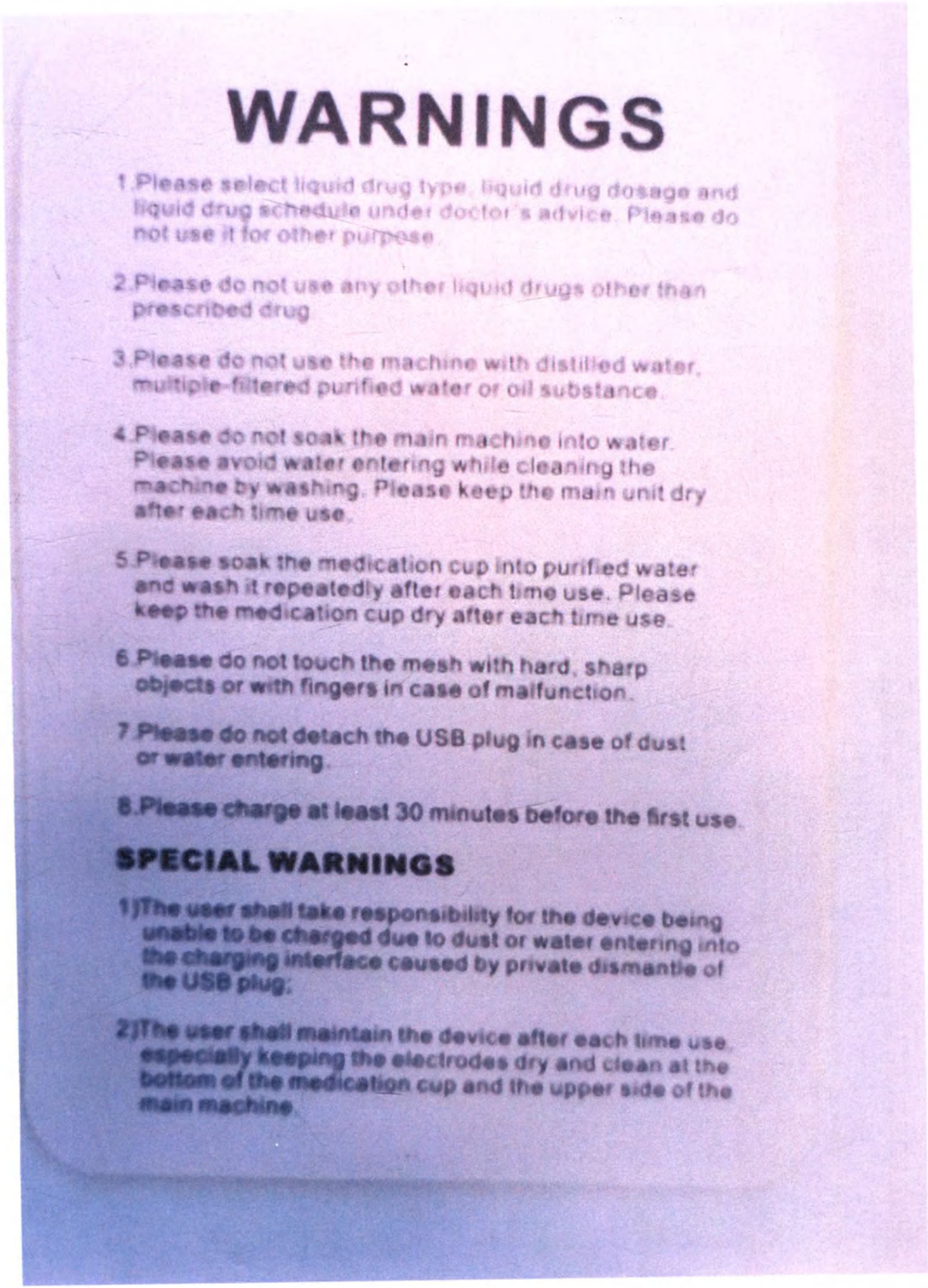
WARNINGS

1.Please select liquid drug type, liquid drug dosage and liquid drug schedule under doctor's advice. Please do not use it for other purpose.

2.Please do not use any other liquid drugs other than prescribed drug.

3.Please do not use the machine with distilled water, multiple-filtered purified water or oil substance.

4.Please do not soak the main machine into water. Please avoid water entering while cleaning the machine by washing. Please keep the main unit dry after each time use.

5.Please soak the medication cup into purified water and wash it repeatedly after each time use. Please keep the medication cup dry after each time use.

6.Please do not touch the mesh with hard, sharp objects or with fingers in case of malfunction.

7.Please do not detach the USB plug in case of dust or water entering.

8.Please charge at least 30 minutes before the first use.

SPECIAL WARNINGS

1)The user shall take responsibility for the device being unable to be charged due to dust or water entering into the charging interface caused by private dismantle of the USB plug;

2)The user shall maintain the device after each time use, especially keeping the electrodes dry and clean at the bottom of the medication cup and the upper side of the main machine.

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